


humantouch.com

800.355.2762





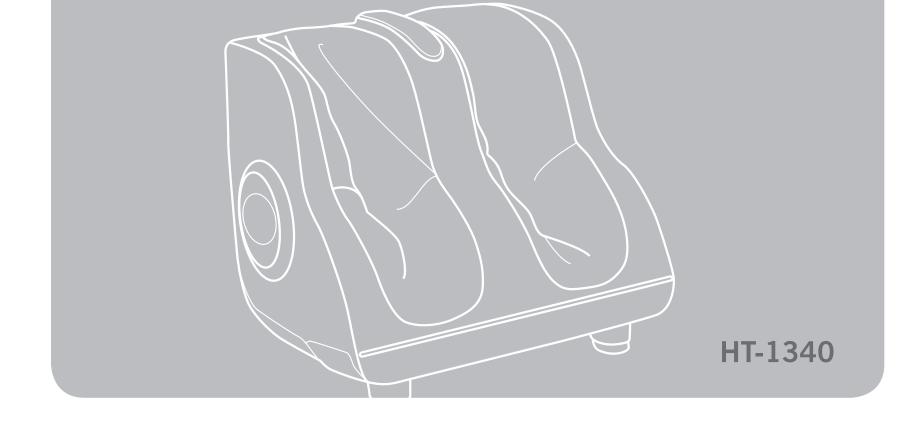

use & care manual

©2009 Human Touch®, LLC. Patents pending and intellectual property rights enforced. No medical claims are warranted or implied by the use of this product. Human Touch is a trademark of Human Touch®, LLC. HT-1340-NA-A0

Human Touch® Partners in Your Well-Being

For decades, Human Touch has been a leader in the field of massage products. We create state-of-the-art products that embrace the latest technology while applying the ancient healing art of massage. Balancing your physical and mental needs requires the right mix of proper nutrition, exercise, and massage. Human Touch massage products are vital components for achieving balance between physical and mental well-being.

To learn more about our company and our products, we invite you to visit us on the Internet at: www.humantouch.com or call our Customer Service Representatives at 800-355-2762.







welcome to human touch®

Your Human Touch® CirQlation™ power foot and calf massager is the perfect bridge between state-of-the-art technology and the ancient healing art of massage. It's a key ingredient for a balanced life—and it's right in your own home!

Your CirQlation™ power foot and calf massager incorporates patented Figure-Eight™ Technology, emulating the same techniques used by massage and chiropractic professionals. It's hard to believe it's not human!

humantouch.com

To ensure that you don't overlook any of your Foot & Calf Massager's features, and that you use it safely, read this manual carefully. Then sit back, relax, and enjoy!

Please send any questions or comments regarding this manual to documentation@humantouch.com

keep it safe

IMPORTANT SAFETY INSTRUCTIONS

When using an electrical appliance, basic precautions should always be followed, including the following:

Read all instructions before using this appliance.



DANGER: To reduce the risk of electrical shock

Always unplug this product immediately after using and before cleaning.



WARNING: To reduce the risk of burns, fire, electric shock or injury to persons:

- Do not put hands or fingers near the massage mechanism or rollers while this product is running. Touching the mechanism can result in injury from the rollers squeezing together.
- If the supply cord is damaged, it must be replaced by the manufacturer or its service agent or a similarly qualified person in order to avoid a hazard
- An appliance should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts.
- Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- Close supervision is necessary when this product is used by, on, or near children, invalids, or disabled persons.
- Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return this product to a service center for examination and repair.
- Do not carry this appliance by supply cord or use cord as a handle.
- Always keep the power cord away from heated surfaces.

- Never operate this product with the air openings blocked. Keep the air openings free of lint, hair, and the like.
- Never drop or insert any object into any opening.
- Do not use outdoors.
- Do not operate where aerosol (spray) products are being used or where oxygen
- To disconnect, turn all controls to the OFF position, then remove plug from outlet.
- Do not massage any area of the body that is swollen, inflamed or covered with skin eruptions.
- Do not use on unexplained calf pain.
- Do not use this product when standing on a damp floor, or when any part of the body is in contact with plumbing or any similar ground.
- Do not stand on or in appliance. Use only while seated.
- Anyone with a medical condition, including those who are or may be pregnant, should consult a physician before using any massage product.

Safety Precautions

- Comfortable clothing is recommended while using this product.
- Do not sit on the controller.
- Do not pull on the AC cord.
- When using this product after it has been in storage, check to see that the product is operating properly prior to use.

SAVE THESE INSTRUCTIONS

Do not attempt to repair this product yourself.

- This product is designed for personal in-home use. Commercial use will void warranty.
- To reduce the risk of electric shock, this appliance has a polarized plug (one blade is wide than the other). This plug will fit in a polarized outlet only one way. If the plug does not fit fully in the outlet, reverse the plug. If it still does not fit, contact a qualified electrician to install the proper outlet. Do not change the plug in any way.

Care and Maintenance

The product is specially designed and constructed to be maintenance free. No periodic lubrication or servicing is required. Your product should only be serviced by an authorized service center. Foreign objects should not be inserted between the massage nodes or the unit housing. Care should be taken not to overload the unit.

STORAGE:

Your product has been carefully designed to give you many years of trouble-free service. Always unplug the unit when not in use. Additionally, if not in use for a long period of time, we recommend that the power cord be coiled and that the unit be placed in an environment free of dust and moisture. Do not store near heat or open flame. Do not leave this product exposed to direct sunlight for long periods of time. This may cause fading or damage.



WARNING: This appliance is equipped with a thermal resetting protection device. This is a safety feature to protect the appliance from overheating. If the appliance should suddenly stop and will not start, turn off the main switch and pull the plug off the socket outlet. Do not operate the appliance for at least 30 minutes. Failure to turn the appliance off may result in the appliance starting unexpectedly when the device cools.

FABRIC CARE:


To maintain your Human Touch product in optimum condition, dust it frequently and/or vacuum it with a soft brush vacuum extension. Faux Suede products can also be brushed with a suede brush. Wipe with a slighty damp cloth as necessary. For tougher stains on non-leather products, use a mild detergent mixed with water. Leather products should be wiped down occasionally with a leather furniture cleaning solution.



WARNING:

Be sure to test any cleaning solution on a small, inconspicuous part of the product to ensure that the color does not fade prior to using the solution on the entire product surface.

get ready

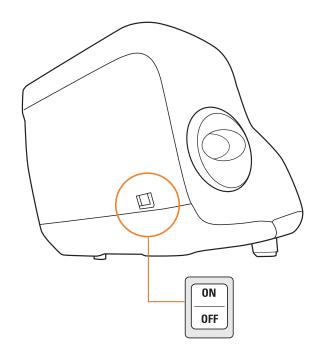



1. Plug it in

Plug the power cord into an appropriate power source.


2. Power on

Turn the master power switch (on the rear lower right-hand corner of the massager) to On.

3. Position your CirQlation power foot and calf massager in front of your favorite chair.



5. Place your feet all the way into the foot and calf massage wells.

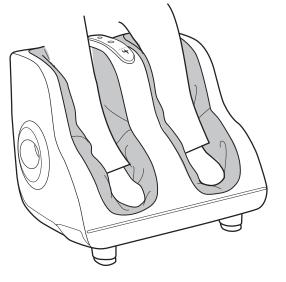


Your CirQlation power foot and calf massager comes fully assembled. Simply plug it in, take a seat and you're set to go!



NOTE: Follow the below steps to remove your feet from the unit in the event of a power failure;


- 1. Lift your knees upward (Fig. 1).
- 2. Rotate your feet sideways and ease them out (Fig. 2).

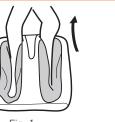




Sitting or standing on your CirQlation power foot and calf massager could result in personal injury and or permanent damage to the massager.

take control

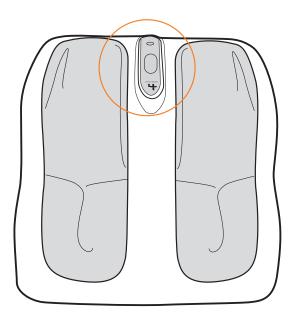

The CirQlation power foot and calf massager automatically limits your massage to an ideal 15 minute interval. To restart, simply re-select your favorite massage setting.

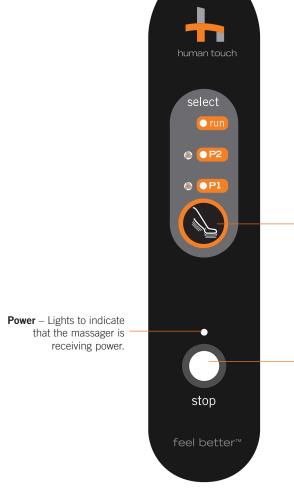

TWO AUTOMATIC PROGRAMS

P1: Gentle Relief

A pulsating kneading massage that improves circulation.

P2: Rapid Refresh

A brisk, effective kneading massage that provides continuous relief.

+ SINGLE SPEED

A continuous cycle of kneading and under-foot massage to relax and regenerate.

Massage Select – Press this button repeatedly to toggle through the fixed speed (run) and auto massage programs (P1, P2) ⊚ . The LED lights to indicate which massage is currently selected.

Stop – Immediately cancels the massage and stops all motion.

Do not massage any single area for more than 20 minutes. If your feet, calves or legs feel uncomfortable, discontinue use immediately.

Your CirQlation power foot and calf massager is very flexible. You can choose a single-speed massage, or select from two auto massage programs, which feature varying speeds.



Patented Figure-Eight™ Technology



Massages calves in an upward, rolling or wave-like" fashion, forcing blood away from the feet and toward the body's core to improve blood flow and wellness.

Non-Skid Floor Protectors -Floor protectors keep unit in place and ensure quiet operation.

Motor-Driven Rollers - Powerful rollers

just under the feet provide a robust

massage to the sole and heals.

limited warranty

Please read the following warranty information for important disclosures before using your new CirQlation power foot and calf massager.

For Targeted Relief Personal Wellness products, the Consumer is fully covered for all labor, parts and shipping costs for the first 90 days of ownership. From day 91 through the remaining Parts Coverage period, the Consumer is responsible for a Repair charge in addition to shipping and handling costs associated with repair or replacement of their unit.

Human Touch CirQlation power foot and calf massager	
Labor, Parts and Shipping	Parts Coverage
First 90 days	1 Year

In the unlikely event that your Human Touch product malfunctions or requires repair, please contact your retailer. For USA customers only, please contact Human Touch Customer Service at (800) 355-2762 between the hours of 6 AM and 5 PM PST or cs@humantouch.com for assistance. You can also send a fax to customer service at (800) 840-2940.

All warranty service requires approval and authorization by Human Touch.

PRODUCT REPLACEMENT WARRANTY PLANS:

Human Touch offers supplemental Product Replacement coverage plans to enhance ownership of your Human Touch product. Replacement Plans are available for purchase directly from Human Touch Customer Service Department. Please call for plan details and pricing.

General Limitations:

- The warranty is only valid within the United States and Canada provided the product has been operated according to the instructions accompanying it and has not been misused or damaged in any way that is beyond the control of Human Touch.
- All service covered by this warranty must be approved by Human Touch and repairs performed by authorized technicians.
- The warranty applies only to the massage chair; accessories or enhancements are not covered.
- Important exclusions apply to cosmetic items such as fabrics, woods, foams, padding, plastics covers, etc.
- Human Touch is not responsible in any way for losses, damages, or inconvenience caused by equipment failure.
- For a complete listing of warranty coverage, exclusions, and limitations, visit humantouch.com.
- Warranties do not apply to rental, business, commercial, institutional, or other non-residential users.
- Proof of purchase (original receipt) is required for all warranty repairs.

Specifications*

oltage: AC110-120V 60 Hz

Power Consumption: 50 Watts Product Weight: 21 lbs.

*Specifications are subject to change without prior notice.

humantouch com