**i.comfort™ Massage Cushion
with Heat**

Model: IBM80010

User Manual

i.comfort™ Massage Cushion with Heat

Model: IBM80010

USER MANUAL

CONTENTS

- Introduction1
- Important Safety Instructions1
- Precautions1
- Overview
 - Massager2
 - Control Panel.....2
- Features2
- Getting Started
 - Turn Massager on2
 - Start Massage.....3
 - Customise Massage..... 3
- Suggestions for Use3
- Maintenance
 - Storing3
 - Cleaning3
- Specifications3
- About Oregon Scientific4
- EU-Declaration Of Conformity4

INTRODUCTION

Thank you very much for purchasing “IBM80010 i.comfort™ Massage Cushion with Heat”. This instruction manual contains useful information on appropriate use and maintenance of this product. In order to benefit from its features, read this manual before first use and save it for future reference.

IMPORTANT SAFETY INSTRUCTIONS

WHEN USING AN ELECTRICAL PRODUCT, ESPECIALLY WHEN CHILDREN ARE PRESENT, ALWAYS FOLLOW BASIC PRECAUTIONS, INCLUDING THE FOLLOWING INSTRUCTIONS: Read all instructions before using these units.

DANGER To reduce the risk of electric shock:

- Unplug the unit from the electrical outlet immediately after using and before cleaning.
- If the unit has been dropped into water, do not pick up or touch it. Switch off at the mains and unplug immediately. Do not operate if it is wet.
- Do not use while bathing or in shower. Keep the unit always dry. Store the unit in a dry and safe place. Keep away from risks of dropping into a sink or a tub. Do not submerge in water.
- Do not insert pins or other metallic fasteners into the unit.
- Keep dry – Do not operate in a wet or moist condition.
- Do not sit or stand on the massage unit.
- To reduce the risk of burns, fire, electric shock, or injury to persons:
 - Keep out of reach of children.
 - Do not leave a plugged unit unattended.
 - Unplug from outlet when not in use and before putting on or taking off parts.
- These units are not intended for use by persons with certain disabilities or lack of knowledge.
- Do not use if you have poor blood circulation. If in doubt, consult with your doctor.
- Do not cover or use the unit under a blanket when operating. Covering it during use may excessively heat the unit which may cause fire, electrocution or injury to persons.
- Stop using the unit when you are feeling drowsy. Do not use when you are sleeping.
- Stop using the unit if it is damaged, wet, tempered, forced or not functioning properly. If in doubt, contact our Helpline.
- Keep the power cord away from heat source.
- Do not lift this unit by pulling the power cord or use the cord as handle.
- Do not drop or insert any object into any opening.
- Do not operate where aerosol (spray) units are being used or where oxygen is being administered.
- Before unplugging, turn the unit off by pressing

the on / off button on the remote control.

- Do not put heavy objects on the unit. Avoid risk of crushing.
- This unit is intended for indoor use only.
- This unit is designed for non-professional use only.
- Use this unit only for its intended use as described in this manual. Do not use provided attachments not recommended by Oregon Scientific.
- Do not drop or insert any object into the unit.
- Save the user manual for future reference.

PRECAUTIONS


Please read all instructions carefully before operating.

- This unit is not intended for medical use or to substitute medical care. It is intended only to provide a luxurious massage.
- Individuals, who are pregnant, fitted with pacemaker, suffering from diabetes, phlebitis and / or thrombosis, or are at an increased risk of developing blood clots due to recent surgery or prolonged bed rest should consult their physician before using the unit.
- Not recommended for use by Diabetics.
- This unit should not be used by any individual suffering from any physical ailment that would limit the user's capacity to operate the controls or who has sensory deficiencies in the lower half of their body.
- Do not leave the unit unattended when children are present.
- This unit is not suitable for children's use. Keep out of reach of children.
- Do not cover the unit with anything when it is in operation.
- Do not obstruct with objects or thrust any part of your body in the moving massage mechanism.
- Do not use the unit for more than 15 minutes at a time. Excessive use may overheat the unit. Should this occur, stop using the unit and wait until it cooled down before the next use.
- Do not use the unit directly on open wounds, swollen or inflamed areas and skin eruptions.
- Do not use this unit before going to bed. The massage has a stimulating effect which may delay sleep.
- Do not use this unit in bed.
- Not recommended for use on timber surfaces as the zipper may damage the wood. Take special care when using on upholstered furniture.
- Do not use in automobiles.
- If you feel ill after use or experience pain in muscles or joints for a long time, discontinue use and consult a physician. Persistent pain could be a symptom of a more serious condition.
- Children should be supervised to ensure that they do not play with the unit.

- Do not use the unit on erupted, swollen, inflamed skin areas, when you have a skin disease, ulcer or varicose veins are present.
- Do not use this unit for more than 15 minutes at a time. Allow 15 minutes rest before next section. Extensive use could lead to the unit's excessive heating and shorter life. Should this occur, discontinue use and allow the unit to cool before operating.
- The unit is equipped with overheating protection circuit, if the power is cut-off, please wait until the unit is cooled down.
- Do not use if you have diabetes, tuberculosis, benign or malignant tumors, phlebitis and thrombosis, or if you have poor circulation. Consult a physician before using if you have any kind of fungal infection, swelling, fracture or persistent pain.
- If you are feeling drowsy, stop using the unit and do not leave the unit operating when you are sleeping.
- Ensure that your hands are dry when operating the controls or removing the plug.
- Do not sit or stand on the massage unit.
- Before unplugging, turn the unit off on the control panel.
- This unit is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the unit by a person responsible for their safety.
- The unit has a hot surface. Persons insensitive to heat must be careful when using the unit.

OVERVIEW

MASSAGER

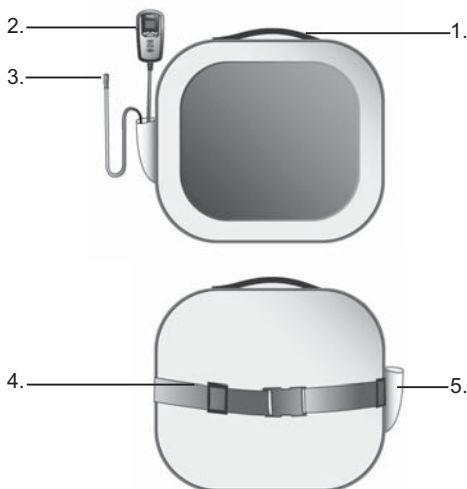1. Carrying strap
2. Remote control
3. Power cord
4. Adjustable strap
5. Remote control storage pouch

REMOTE CONTROL

1. LCD screen: indicate selection / time remaining
2. ↻ / ↺ : Adjust massage direction (clockwise / anticlockwise)
3. **ON / OFF** : Turn i.comfort on / off
4. **+** / **-** : Adjust massage intensity (3 programmable levels)
5. **HEAT**: Turn heating function on / off

FEATURES

- Optional heat massage function with light indication
- Versatile massager relieves muscle fatigue where you need it most, from your upper back and shoulders, thighs, calves, lower back and feet
- Customise your 4 point kneading massage with a touch of a button (3 intensities and 2 motions to choose from)
- Bright contemporary design with sleek ergonomic remote control for urban living
- Compact design with comfortable hand strap for easy portability
- Lightweight with adjustable strape for easy attachment to most seats - Just plug it in and relax
- Easy to clean synthetic leather finish
- 15 minute energy saving timer with auto shut off

GETTING STARTED

TURN MASSAGER ON

Plug adapter into socket and connect to power cord. 'OF' on LCD Screen indicates massager is in standby mode.

NOTE

- Make sure the adapter is not obstructed and is easily accessible to the personal massager.
- To completely disconnect from power, the adapter should be disconnected from the power cord by pulling it out while holding the plug. Do not disconnect with wet hands or yank cord out.
- The massager and adapter should not be exposed to wet conditions. No objects filled with liquid, such as vases, should be placed on the personal massager and its adapter.

START MASSAGE

1. Place the massager behind your back, under your feet or legs.

NOTE

Adjustable straps provided can be used to affix the massager to seat.

2. Press **ON/OFF** button to start operating the massager. '15' indicates massage has begun and there are 15 minutes remaining.


CUSTOMISE MESSAGE

1. Press  /  to select desired massage direction ('C1' indicates clockwise, 'C2' indicates anticlockwise).
2. Press **+** / **-** to select desired massage intensity ('P1' indicates normal, 'P2' indicates medium, 'P3' indicates strong).
3. Press **HEAT** button to activate the heat function. To turn the heat off, Press the **HEAT** button again. We advice not to exceed 15-minute massage duration at a time.


NOTE

The unit will be stopped automatically when overheat.

SUGGESTIONS FOR USE

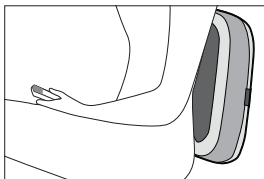
Lower Back

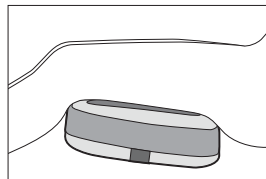
Thighs

Shoulder

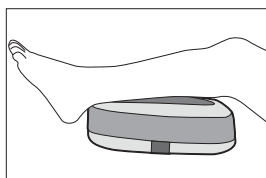Calves

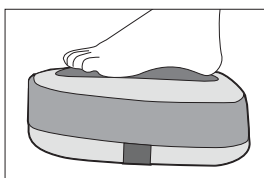Feet

NOTE

This product is not intended to be used as a pillow or headrest. Do not place around head area and / or use as head massager.

WARNING

This unit can only withstand a maximum of 30 kg pressure. Do not sit or rest heavy objects on unit. For safety reasons, unit will automatically turn off if over 30 kg pressure is exerted on it.

MAINTENANCE

STORING

Place the unit in its box or in a safe, dry, cool place. Refrain from sharp edges or pointed objects that might cut or puncture the fabric surface. To avoid breakage, do not wrap the power cord around the unit. Do not hang the unit by the power cord or the cord of the remote control unit.

For your travelling comfort and ease of storage, a woven bag has been provided.

CLEANING

Unplug the unit and allow it to cool down before cleaning. Clean only with a soft, slightly damp sponge. Never allows water or any other liquids to come into contact with the unit.

- Do not immerse in any liquid.
- Never use abrasive cleaners, brushes, gasoline, kerosene, glass / furniture polish or paint thinner to clean.
- Do not remove or unzip the unit cover to clean the internal parts.

Do not attempt to repair the i.comfort™ Massage Cushion with Heat. There are no user serviceable parts.

SPECIFICATIONS

| TYPE | DESCRIPTION |
|-------------------|--|
| MAIN UNIT | |
| L x W x H | 410 (L) x 350 (W) x 135 (H) mm (30.7in x 15.7in x 3.54in) |
| Weight | 2.9 Kg |
| Power Consumption | 36W |
| Voltage | Input 100-240V ~ 50 / 60Hz Output 24V / 1500mA |

NOTE

The technical specifications for this product and the contents of the user manual are subject to change without notice.

ABOUT OREGON SCIENTIFIC

Visit our website (www.oregonscientific.com) to learn more about Oregon Scientific products. If you're in the US and would like to contact our Customer Care department directly, please visit:

www2.oregonscientific.com/service/support

OR

Call 1-800-853-8883.

For international inquiries, please visit:

www2.oregonscientific.com/about/international

EU-DECLARATION OF CONFORMITY

Hereby, Oregon Scientific, declares that this i.comfort™ Massage Cushion with Heat (model: IBM80010) is in compliance with the essential requirements and other relevant provisions of Directive 2004/108/EC. A copy of the signed and dated Declaration of Conformity is available on request via our Oregon Scientific Customer Service.

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