


Flex Targeted Massager

SPOT-ON RELIEF FOR STIFF, ACHING MUSCLES

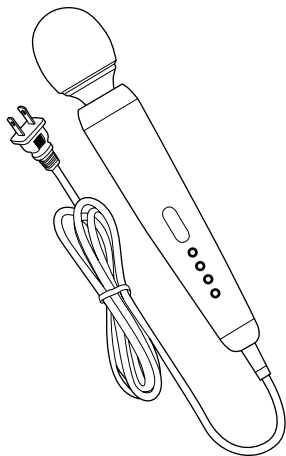

CAUTION

- TO REDUCE THE RISK OF ELECTRIC SHOCK, DO NOT REMOVE COVER. THERE ARE NO SERVICEABLE PARTS INSIDE.
- TO REDUCE THE RISK OF FIRE OR ELECTRIC SHOCK, DO NOT EXPOSE THIS UNIT TO RAIN OR MOISTURE.



The lightning flash with arrow-head symbol within an equilateral triangle is intended to alert the user to the presence of uninsulated "dangerous voltage" within the unit's enclosure that may be of sufficient magnitude to constitute a risk of electric shock.



The exclamation point within an equilateral triangle is intended to alert the user to the presence of important operating and maintenance (servicing) instructions in the literature accompanying the unit.

IMPORTANT SAFETY INSTRUCTIONS

All of the safety and operating instructions should be read, adhered to and followed before the unit is operated.

SAVE THESE INSTRUCTIONS

DANGER

To reduce the risk of electric shock, burns, fire or injury:

1. Always unplug the product from the electrical outlet immediately after using and before cleaning.
2. Do not use while bathing or in a shower.
3. Do not place or store product where it can fall or be pulled into a tub or sink.
4. Do not place in or drop into water or other liquid.
5. Do not reach for product that has fallen into water. Unplug it immediately.

WARNING

1. Do not leave product unattended when plugged in. Unplug from outlet when not in use.
2. This product should not be used on or by children.
3. Do not use attachments not recommended by manufacturer.
4. Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water.
5. Keep cord away from heated surfaces.

6. Never block the air openings of the product or place it on a surface, such as a bed or couch, where the air openings may be blocked. Keep the air openings free of lint and hair.
7. Never use while sleeping.
8. Do not use outdoors.
9. Never drop or insert objects into any openings.
10. Do not operate where aerosol products are being used or where oxygen is being administered.
11. Do not use on or near eyes or other highly sensitive areas.
12. Do not operate under a blanket or pillow. Excessive heating can occur and cause fire, electric shock or injury.
13. Do not carry this product by its power cord or use the cord as a handle.
14. Do not use directly on one concentrated area of the body. It is meant to be used in a continuous motion, constantly changing the massage area.
15. To disconnect, turn all controls to the off position, and then remove the plug from the outlet.
16. Use this appliance only for its intended use as described in this manual. No serviceable parts inside.


17. This product is designed for external body use only.
18. For adult use only.

CAUTION

- Consult your doctor before using this product.
- Individuals with pacemakers should not use this product.
- Never leave the product operating or unattended, especially if children are present.
- Never use the product on open skin eruptions or on swollen or inflamed areas.
- Do not operate after taking painkillers, sedatives or alcohol. Consult your doctor if under these medications.
- Do not operate for more than 20 minutes continuously. It is recommended to turn off the product and allow it to rest (cool) for 40 minutes before re-using. This will prolong the life of the product.
- The product is designed as a non-professional massager to soothe tired and aching muscles. It is not intended to be a substitute for medical attention. If your ailments and symptoms persist or worsen after use, contact your doctor immediately.

LOCATION OF PARTS AND CONTROLS

1. Vibrating Head
2. Main Body
3. Power/Speed Button
4. LED Speed Indicator (4 Speeds)
5. Power Cord

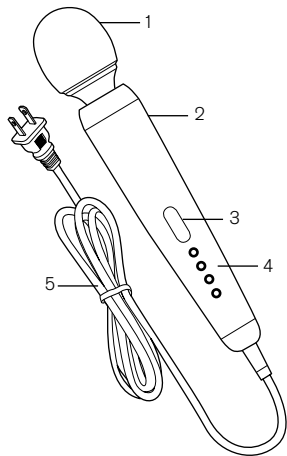

OPERATION

This massager is equipped with a single button for easy operation.

1. Press the Power/Speed Button to turn the massager ON to the 1st level.
2. Press the Power/Speed Button to increase speeds through the four intensity levels. The LED SPEED INDICATOR lights will light up to let you know what level you are using.

- 1 light - level 1
- 2 lights - level 2
- 3 lights - level 3
- 4 lights - level 4

3. Press the Power/Speed Button to turn OFF the massager from level 4.
4. Gently apply to desired area for targeted relief.

Note: Massager automatically shuts off after 20 minutes of continuous use.

IMPORTANT: To reduce the risk of electric shock, this product has a polarized plug (one blade is wider than the other). This plug will fit in a polarized outlet only one way. If the plug does not fit fully in the outlet, reverse the plug. If it still does not fit, contact a qualified electrician to install the proper outlet. **DO NOT CHANGE THE PLUG IN ANY WAY.**

CARE AND MAINTENANCE

1. Unplug the massager from the electrical outlet before cleaning.
2. To prevent germs from spreading, carefully clean before and after use using a soft, damp cloth.
3. DO NOT clean this massager with alcohol or any liquid chemical cleansers.
4. Ensure the massager is completely dry before use.
5. Store in a cool, dry place.
6. Do not twist cord around product while storing. Gather cord together and secure with included hook and loop strap.

SPECIFICATIONS

Voltage: 120 V, 60Hz

Node Vibration (approximately):

Level 1: 3,300 rpm

Level 2: 4,200 rpm

Level 3: 5,000 rpm

Level 4: 6,000 rpm

Weight: 1.5 lbs (approximately)

Dimensions: 13.40" x 1.25" x 2.38"

Cord Length: 7 ft.

ONE (1) YEAR LIMITED WARRANTY

Brookstone® warrants this product against defects in materials and/or workmanship under normal use for a period of ONE (1) YEAR from the date of purchase by the original purchaser ("Warranty Period"). If a defect arises and a valid claim is received within the Warranty Period, at its option, Brookstone will either 1) repair the defect at no charge, using new or refurbished replacement parts, or 2) replace the product with a new product that is at least functionally equivalent to the original product, or 3) provide a store credit in the amount of the purchase price of the original product. A replacement product or part, including a user-installable part installed in accordance with instructions provided by Brookstone, assumes the remaining warranty of the original product. When a product or part is exchanged, any replacement item becomes your property and the replaced item becomes Brookstone's property. When a store credit is given, the original product must be returned to Brookstone and becomes Brookstone's property.

Obtaining Service: To obtain warranty service, call Brookstone Limited Warranty Service at 1-800-292-9819. Please be prepared to describe the product that needs service and the nature of the problem. A purchase receipt is required. All repairs and replacements must be authorized in advance. Service options, parts availability and response times will vary. You are responsible for delivery and the cost of delivery of the product or any parts to the authorized service center for replacement, per our instructions.

Limits and Exclusions: Coverage under this Limited Warranty is limited to the United States of America, including the District of Columbia and the U.S. Territories of Guam, Puerto Rico, and the U.S. Virgin Islands. This Limited Warranty applies only to products manufactured for Brookstone that can be identified by the "Brookstone" trademark, trade name, or logo affixed to them or their packaging. The Limited Warranty does not apply to any non-Brookstone products. Manufacturers or suppliers

other than Brookstone may provide their own warranties to the purchaser, but Brookstone, in so far as permitted by law, provides these products "as is." This warranty does not apply to: a) damage caused by failure to follow instructions relating to product's use or the installation of components; b) damage caused by accident, abuse, misuse, fire, floods, earthquake or other external causes; c) damage caused by service performed by anyone who is not a representative of Brookstone; d) accessories used in conjunction with a covered product; e) a product or part that has been modified to alter functionality or capability; f) items intended to be periodically replaced by the purchaser during the normal life of the product including, without limitation, batteries or light bulbs; g) any product sold "as is" including, without limitation, floor demonstration models and refurbished items; or h) a product that is used commercially or for a commercial purpose.

BROOKSTONE SHALL NOT BE LIABLE FOR INCIDENTAL OR CONSEQUENTIAL DAMAGES RESULTING FROM THE USE OF THIS PRODUCT, OR ARISING OUT OF ANY BREACH OF THIS WARRANTY. TO THE EXTENT PERMITTED BY APPLICABLE LAW, BROOKSTONE DISCLAIMS ANY AND ALL STATUTORY OR IMPLIED WARRANTIES, INCLUDING, WITHOUT LIMITATION, WARRANTIES OF MERCHANTABILITY, FITNESS FOR A PARTICULAR PURPOSE AND WARRANTIES AGAINST HIDDEN OR LATENT DEFECTS. IF BROOKSTONE CANNOT LAWFULLY DISCLAIM STATUTORY OR IMPLIED WARRANTIES, THEN TO THE EXTENT PERMITTED BY LAW, ALL SUCH WARRANTIES SHALL BE LIMITED IN DURATION TO THE DURATION OF THIS EXPRESS WARRANTY.

Some states disallow the exclusion or limitation of incidental or consequential damages or how long an implied warranty lasts, so the above exclusions or limitations may not apply to you. This warranty gives you specific legal rights and you may also have other rights, which vary from state to state.