<b>DK</b>	Dagslys Lampe .....	2
<b>SE</b>	Dagsljuslampa .....	5
<b>NO</b>	Dagslys Lampe .....	8
<b>FI</b>	Päivänvalolamppu .....	11
<b>UK</b>	Daylight lamp .....	14
<b>DE</b>	Tageslichtlampe .....	17
<b>PL</b>	Lampa dzienna .....	21

## INTRODUKTION

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## INTRODUCTION

To get the best out of your new appliance, please read this user guide carefully before using it for the first time. Take particular note of the safety precautions. We also recommend that you keep the instructions for future reference, so that you can remind yourself of the functions of your appliance.

## SAFETY MEASURES

### General

- Incorrect use of this appliance may cause personal injury and damage.
- Use for its intended purpose only. The manufacturer is not responsible for any injury or damage resulting from incorrect use or handling (see also Guarantee Terms).
- For domestic use only. Not for outdoor or commercial use.
- For connection to 230 V, 50 Hz only.
- Do not immerse the appliance, cord or plug in water and make sure no water enters the appliance.
- Do not use in humid surroundings.
- Never touch the appliance, cord or plug with wet or damp hands.
- Do not expose to direct sunlight, high temperatures, humidity, noise or corrosive substances.
- Keep the appliance under constant supervision while in use. Keep an eye on children.
- Do not cover.
- Use under temperature 0-35 °C.
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.

### Cord and plug

- Keep the cord and appliance away from heat sources, hot objects and naked flames.
- Check that it is not possible to pull or trip over the cord or any extension cord.
- Switch off the appliance, remove the plug from the socket when the appliance is not in use and before cleaning or replacement of the fluorescent tubes. Avoid pulling the cord when removing the plug from the socket. Instead, hold the plug.
- Handle fluorescent tubes carefully as they may splinter and cause injuries to your hands.
- Ensure that the cord is fully extended.
- The cord must not be twisted or wound around the appliance.
- Check regularly for damage to appliance, cord and plug. Do not use if any parts are damaged or if the appliance has been dropped or damaged in any way whatsoever.


- If the appliance, cord or plug have been damaged, have them inspected and if necessary repaired by an authorised repair engineer. Never try to repair the appliance yourself. Please contact the store where you bought the appliance for repairs under guarantee. Unauthorised repairs or modifications to the appliance will invalidate the guarantee.

## Using lightbaths

- Do not use in rooms without some other light source. For use as a light supplement only, not illumination.
- Do not look directly at the appliance for long periods. Follow the user guide below.
- Initial use may cause tired eyes and/or headache. This is normal, and such effects will soon disappear. If symptoms persist, consult your doctor.
- Consult your doctor before using this appliance if you
  - suffer or have suffered from severe depression
  - are hypersensitive to light
  - suffer from any type of eye disease
  - use particular forms of medicine such as psychopharmacia or malaria pills
  - suffer from high blood pressure.
- Consult your doctor, if the use of the appliance results in hyperactivity, exhaustion, depression, restlessness, difficulties in sleeping and/or other adverse effects.

## KEY TO MAIN COMPONENTS OF THE APPLIANCE

1. Top cover
2. Screen
3. Rear part
4. Stand
5. On/Off button
6. Front base
7. Rear base
8. Cord with plug

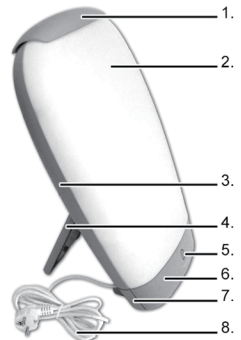Wall bracket for mounting  
(not shown on the illustration)

## USING THE APPARATUS

The appliance is used for lightbaths that counteract winter fatigue, which can cause the following symptoms:

- fatigue
- bad mood or moodiness
- lack of energy
- increased need for sleep
- lack of concentration
- lack of appetite

Over 60% of the population are believed to experience one or more symptoms of winter blues due to the low intensity of the winter sun.

The most effective time of day to take a light bath is between 6 am and 8 pm for a recommended period of 2 hours per day.

You can take lightbaths as often as you like but the treatment is most effective if you take lightbaths of the below-mentioned duration at least five days in a row.

### Positioning the appliance

- Unfold the stand (4) at the back of the appliance and place the appliance on a table.  
or
- Mount the appliance on the wall using the enclosed wall brackets. Check before mounting that the wall and the screws can carry the weight of the appliance.

### Use

- Insert the plug (8) in the mains socket and switch on the power.
- Switch on the on/off button (5).
- Sit down in front of the appliance in a position where you can comfortably see the light and with your face 20-60 cm from the screen (2). You do not need to sit directly in front of the screen.
- The duration of the lightbath depends on your distance from the screen:

Distance to screen	Max. length of light bath
60 cm	Approx. 2 hours
40 cm	Approx. 1 hour
20 cm	Approx. ½ an hour

- The lightbath can be interrupted and you can do other things at the same time (e.g. read or watch television) as long as you sit within 60 cm from the appliance.
- Do not look directly into the light constantly, but look at it regularly.
- Always turn off the appliance using the on/off button, and unplug it, when it is not in use.
- Allow the appliance to cool for at least 15 minutes before putting it away or moving it.

TIP: Carry the appliance by the carrying handle on the top of the rear part (3).


### REPLACEMENT OF FLUORESCENT TUBES AND STARTING SWITCH

The lifetime of the fluorescent tubes is approx. 8,000 hours of normal use. If a tube fails, it must be replaced.

- **IMPORTANT!** Discard used fluorescent tubes in accordance with the rules for environmentally harmful waste, e.g. by taking them to a waste recovery site.

### How to dismantle the appliance

- **IMPORTANT!** Always make sure to switch off the appliance and remove the plug from the mains socket before taking it apart to avoid the risk of electric shocks.
- Place the appliance on a non-scratching surface, screen down.
- Remove the top cover (1) by unfastening it with your fingers, first in one side and then in the other.
- Remove the two screws at the top of the stand using a Phillips screwdriver.
- Remove the two screws in the rear base (7), and carefully loosen the rear and front parts of the base (6) from the rest of the appliance.
  - **IMPORTANT!** Very carefully detach the two base parts, since they are connected to the rest of the appliance with cords that must not be damaged.
- Remove the 9 screws along the side of the rear part.
  - The screws of the appliance may be of different types. Therefore, pay attention to which screws belong to the different screw holes in order to place them correctly when reassembling the appliance.
- Carefully remove the rear part, turn it upside down and put it down in front of you.
- Replace the starting switch (if applicable):
  - Turn the starting switch (b) anti-clockwise and remove it.
  - Insert a new starting switch of the same type and turn it clockwise until it locks into place.
- Replace the fluorescent tubes:
  - Remove the locking clip (d) by pushing the ends upwards until they click into place and remove the clip.
  - Grip the socket of the fluorescent tube (a) and pull the tube out of the socket (towards the top of the appliance).
  - Insert a new fluorescent tube in the socket as far as it goes until it clicks into place.
  - Press the locking clip onto the holder (c) until it clicks into place. Check that the fluorescent tube is fitted correctly.
  - **IMPORTANT!** Replace the fluorescent tubes with tubes of the same type. Always replace both tubes at the same time.
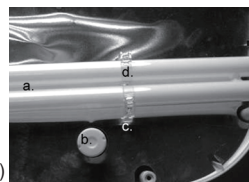
### How to reassemble the appliance

- Place the rear part on the screen, so that the screw holes along the side of the rear part align with the holes on the screen. The two parts click together.
- Assemble the front base (with the on/off button) and make sure the cords are not caught between the parts.
- Attach the rear base and make sure the stand is correctly aligned with the screw holes at the back of the appliance.
- Screw the two parts of the base together and fasten the stand with the two screws at the top.

- Place the rest of the screws in the screw holes along the side of the appliance and fasten them.
- Assemble the top cover.

## CLEANING

When cleaning the appliance, you should pay attention to the following points:

- Remove plug from the mains socket, and allow the appliance to cool down before cleaning it.
- Never submerge the apparatus in water, and no water must be allowed to get into it.

### External cleaning

- Clean the appliance by wiping it with a damp cloth. A little detergent can be added if the appliance is heavily soiled.
- Never use a scouring pad, steel wool or any form of strong solvents or abrasive cleaning agents to clean the appliance, as those may damage the outside surfaces of the appliance.
- IMPORTANT! Allow to dry fully after cleaning before using again. Do not use appliance if damp.

### Internal cleaning

- Remove the screen as described above, wipe off reflector, inside of screen and tubes with a dry, soft cloth. Do not use any kind of cleaning agents.

## SPECIFICATIONS

Model	WL1	WL3
Voltage	220-240 V	220-240 V
Wattage	72 W	72 W
Hz	50	50
LUX	10.000	10.000
Lamp type	FSL-YDW36-HRR	Philips-PL-L36W/865/4P

## INFORMATION ON DISPOSAL AND RECYCLING OF THIS PRODUCT

Please note that this Adexi product is marked with this

symbol:

This means that this product must not be disposed of along with ordinary household waste, as electrical and electronic waste must be disposed of separately.

In accordance with the WEEE directive, every member state must ensure correct collection, recovery, handling and recycling of electrical and electronic waste. Private households in the EU can take used equipment to special recycling stations free of charge. In certain member states, used apparatus can be returned to the dealer where they were bought on the condition you buy new products. Contact your retailer, distributor or the municipal authorities for further information on what to do with electrical and electronic waste.

## GUARANTEE TERMS

The guarantee does not apply:

- if the above instructions are not followed
- if the appliance has been interfered with
- if the appliance has been mishandled, subjected to rough treatment, or has suffered any other form of damage
- if faults have arisen because of faults in your electricity supply.

Due to the constant development of our products in terms of function and design, we reserve the right to make changes to the product without prior warning.

## FREQUENTLY ASKED QUESTIONS

If you have any questions regarding use of the appliance and cannot find the answer in this user guide, please try our website at [www.adexi.eu](http://www.adexi.eu). Go to the "Consumer Service" menu, click on "Question & Answer" to view the most frequently asked questions.

You can also see contact details if you need to contact us for technical issues, repairs, accessories and spare parts.

## IMPORTER

Adexi Group  
[www.adexi.eu](http://www.adexi.eu)

We cannot be held responsible for any printing errors.

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