
TANITA

Monitoring Your Health

Digital Lithium

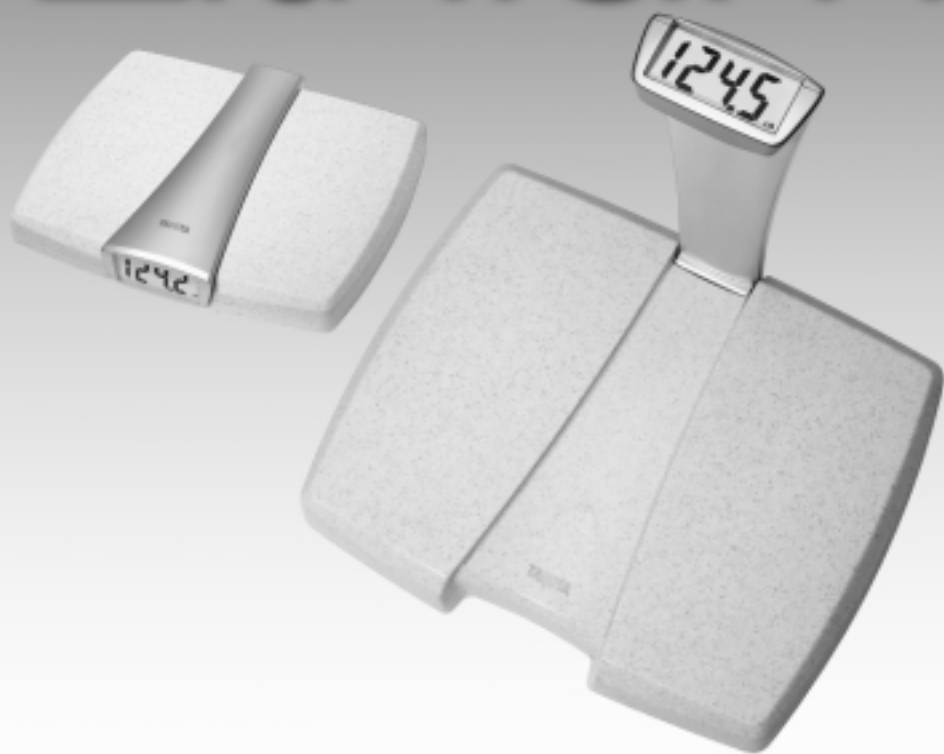
Operating Instructions – HD-335

OPERATING INSTRUCTIONS

Digital Lithium Scale – HD-335

- 1 Measuring Platform
- 2 LCD Display
- 3 Battery Compartment on Bottom of Scale
- 4 Weight Mode Switch on Bottom of Scale
- 5 Display Arm

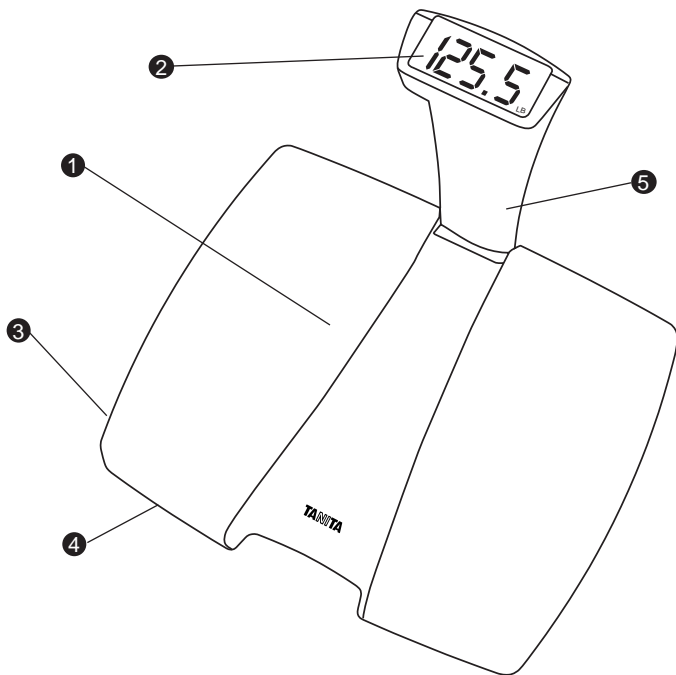Specifications

Maximum Capacity: 350 lb (160 kg)

Graduation: 0.2 lb (0.1 kg)

Power Supply: 1 x 3 volt CR2032 Replaceable Lithium Battery (Included)

UK Specifications

Maximum Capacity: 25 st (160 kg) (350 lb)

Graduation: 1/4 lb (0.1 kg) (0.2 lb)

Power Supply: 1 x 3 volt CR2032 Replaceable Lithium Battery (Included)

Features

- Digital accuracy and reliability.
- Durable construction, non-skid surface.
- “Lo” message appears when battery needs replacing.

Before Use

Select your preferred weight mode using the Weight Mode Switch.

Placing Scale

The scale should be used on a hard even surface for accurate weighing. Not recommended for use on carpeted surfaces. Placing scale on a very soft, deep pile rug may affect weighing accuracy.

Instructions for Use

Press firmly with foot on closed arm to activate display & open arm. If the arm is already up but display does not activate, tap the scale on lower right with your foot to activate display.


Tap the scale with your foot. When "00" shows on the display, step on the scale (Fig. 1).

Stand still while the LCD flashes "00" for 3 seconds. The LCD will display your weight (Fig. 2). If you step off from the scale, the LCD will display your weight for 2 seconds then the scale will automatically shut off. If you remain on the scale, the LCD will display your weight for 10 seconds then automatically shut off.

Note: If you step on the scale without tapping the platform, the display will show "Err" for 2 seconds then automatically shut off.

Once the scale displays "00", if you do not step on the scale in 5 seconds, the scale will show "Err" then automatically shut off.

"EE" will appear on the display when the scale's weight capacity is exceeded.

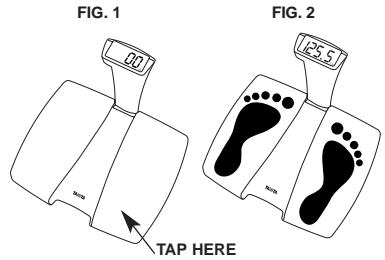
Care and Maintenance

- Place the scale in an area free from heating equipment, high humidity, or extreme temperature changes.
- Never submerge the scale in water or spill chemical liquids on it. To clean, use a dampened cloth and mild detergent.
- Avoid excessive impact or vibration to the unit.

Replacing the Batteries

Lithium batteries can be typically used for many years before replacement is necessary. When the low battery sign, LO, appears on the LCD, please replace with one new lithium battery.

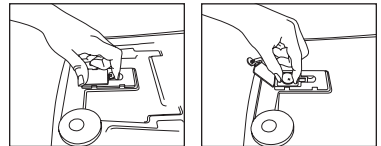Open the battery cover on the back of the measuring platform. Insert the new battery as indicated.

Note: Be sure that the polarity of the batteries is set properly.

If you do not intend to use this unit for a long period of time, it is advisable to remove the batteries before storage.

After Sales Service

Please contact the Tanita office nearest you regarding operation or trouble, request for repair, and/or placement of an order for parts for products.

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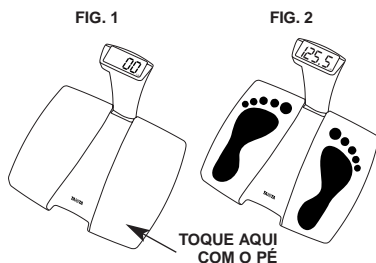
Stand still while the LCD flashes “00” for 3 seconds. The LCD will display your weight (Fig. 2). If you step off from the scale, the LCD will display your weight for 2 seconds then the scale will automatically shut off. If you remain on the scale, the LCD will display your weight for 10 seconds then automatically shut off.

Note: If you step on the scale without tapping the platform, the display will show “Err” for 2 seconds then automatically shut off.

Once the scale displays “00”, if you do not step on the scale in 5 seconds, the scale will show “Err” then automatically shut off.

“EE” will appear on the display when the scale’s weight capacity is exceeded.

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Tanita France S.A.
Tanita Health Equipment H.K. Ltd.
Tanita Europe GmbH
Tanita UK Ltd.