REMINGTON®

SpaTherapy
COLLECTION

Deluxe Paraffin Spa

FEATURES:

- Aromatherapy wax soothes and renews hands, elbows or feet (up to men's size 13)
- Convenient Timer - Paraffin is ready to use when you are
- Variable Heat Control customizes your Paraffin Treatment
- Heel Rest supports the leg and foot while paraffin dries

CONGRATULATIONS...

We at Remington wish to thank you for supporting Remington Products Company, L.L.C. with your purchase. As with all Remington products, your new Paraffin Spa has been manufactured with the utmost care and attention to detail and performance to ensure your complete satisfaction.

Relax and refresh as you escape from everyday stresses and strains in a warm bath of aromatic paraffin wax. Just dip your hand in and out of the Paraffin Spa several times to form a "soft" glove of paraffin wax, which releases therapeutic and moisturizing properties deep into the skin and joints. Immediately you can start to feel the heat actually transferring into your skin. Fifteen brief minutes of indulgence are all it takes to refresh, rejuvenate and invigorate your mind and body.

The Remington Paraffin Spa helps to:

- Soothe and renew hands, feet, elbows and knees
- Moisturize dry, cracked hands, feet, elbows and knees
- Ease stiff joints

IMPORTANT SAFETY INSTRUCTIONS

When using an electrical appliance, especially when children are present, basic safety precautions should always be followed, including the following.

**READ ALL INSTRUCTIONS
BEFORE USING.**

**DANGER
KEEP AWAY FROM WATER**

To reduce the risk of electric shock:

1. Always unplug this appliance from the electrical outlet immediately after using. See cleaning instructions.

2. Do not use while bathing or in the shower.
3. Do not place or store this appliance where it can fall or be pulled into a tub, sink or dropped into water or other liquid.
4. Do not place or drop into water or other liquid.
5. Do not reach for an appliance that has fallen into water, unplug it immediately.

WARNING

To reduce the risk of burns, electric shock, fire, or injury to persons:

1. This product is for household use only.
2. This product should never be left unattended when plugged in. Always unplug from outlet before adding wax and after using.
3. Close supervision is necessary when this product is used by, on, or near children, invalids, or disabled persons.
4. Use this product only for its intended use as described in this manual. Do not use attachments not recommended by Remington.
5. Never operate this product if it has a damaged cord or plug; if it is not working properly; or if it has been dropped or damaged. Return the product to Remington Customer Service for repair or replacement.
6. Keep the cord away from heated surfaces.
7. Never use while sleeping.
8. Never drop or insert any foreign object into the Paraffin Spa chamber.
9. Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
10. Do not use an extension cord with this appliance.
11. This appliance has a polarized plug (one blade is wider than the other). As a safety feature, this plug will fit in a polarized outlet only one way. If the plug does not fit fully in the outlet, reverse the plug. If it still does not fit, contact a qualified electrician. Do not attempt to defeat this safety feature.
12. Do not apply hot wax near eyes or other heat sensitive areas, to any skin injury, eruption, or surface growth.
13. Use this appliance on a level surface.
14. Do not let the cord hang over the edge of the table or counter where it can be pulled on by children or tripped over unintentionally.
15. To disconnect, turn all controls to the off position, then remove plug from the outlet. Allow to cool before handling, cleaning, or storing the appliance.
16. This appliance has no serviceable parts.
17. Do not plug or unplug this product when hands are wet.
18. Do not carry this product by the supply cord or use cord as a handle.
19. Always store this product in a moisture free area.
20. Discontinue use of this appliance if discomfort, irritation, or pain results. Consult your doctor.
21. Do not stand on or in the Paraffin Spa chamber. Use only when seated.
22. Do not pour melted paraffin from the Paraffin Spa chamber.
23. Never use sharp instruments to remove the solidified paraffin from the unit.
24. Do not put water in the Paraffin Spa chamber or immerse the unit in water.
25. Do not allow the level of paraffin in the Paraffin Spa chamber to fall below the "min" line, or above the "max" line. With the unit filled to the "max" line, it takes approximately 8–10 uses (5 dips per foot per use) for the wax to fall below the "min" line.
26. Do not heat paraffin on a range, stove, or in a microwave, or in plastic packaging material.
27. Do not use any additives in the paraffin.
28. **CAUTION: WAX MUST BE IN THE UNIT BEFORE TURNING ON OR DAMAGE MAY OCCUR.**

PLEASE NOTE: This appliance should not be used in the presence of open cuts or wounds, inflammatory skin conditions, growth peripheral vascular disease where circulation is impaired, acute inflammation, or where sensation of the extremity is reduced or absent. If there is any question about peripheral vascular disease or decreased sensation of the extremities, consult a physician or physical therapist before using the Paraffin Spa. Do not use the Paraffin Spa if you have any dermatological conditions, consult a physician first. Discontinue use if dermatitis due to paraffin sensitivity occurs. Discontinue use if wax feels too hot or cool, which could indicate health problems with the user.


This product has no pharmaceutical effect. This paraffin is designed to be used with paraffin heating baths for heat transfer purposes only.

SAVE THESE INSTRUCTIONS

OPERATING INSTRUCTIONS

How to Fill the Spa

1. Remove the lid. Set all documents aside.
2. Place the unit on a hard, flat, heat-resistant surface.
3. Remove all plastic packaging material from the paraffin wax then place 6lbs. of paraffin wax in Paraffin Spa chamber. Do not use any other wax except Remington wax with this product.
4. Make sure unit is "off". Plug unit into an electrical outlet of proper voltage. When the unit is plugged in, the word "off" will show on the display.
5. This product has a polarized plug (one blade is wider than the other). As a safety feature, this plug will fit in a polarized outlet only one way. If the plug does not fully fit in the outlet, reverse the plug. If it still does not fit in the outlet, contact a qualified electrician to install a polarized outlet. Never try to defeat this safety feature.
6. Replace the lid to retain heat.

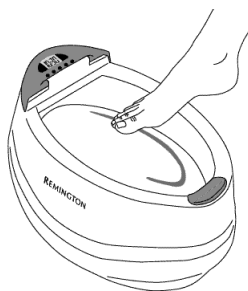How to Use the Spa:

Once you have completed the "How to Fill the Spa" instructions, you can either turn the unit on by pressing the on/off button to commence melting immediately. Or, you can set the timer to commence melting at a later point in the day.

How to Commence Melting Immediately:

1. Turn the unit on by pressing the power "On/Off" button. When the unit is turned on, the temperature will default to medium, and the word "MED" will show on the display.
2. You can switch to the "HIGH" temperature setting by pressing the "TEMP" button once. You can switch to the "LOW" temperature setting by pressing the "TEMP" button twice. Or, you can switch back to the "MED" (default) temperature setting by pressing the "TEMP" button a third time.
3. **It is recommended that the wax be melted on "HIGH" temperature to achieve the fastest melt time (approximately 3 hrs.)**

4. When the wax is visibly melted (approximately 3 hrs. on "HIGH" temperature), the Deluxe Paraffin Spa is ready to use (see "How to dip feet" instructions). Once wax is melted, do not move unit as spillage may occur!
5. Everyone's sensitivity is different. So, if the temperature feels too hot, simply press the "TEMP" button to reach the desired temperature (LOW-MED-HIGH). It will take approximately one hour (from HIGH temperature) for the Paraffin Spa to reach the MED temperature, and two hours (from HIGH temperature) for the Paraffin Spa to reach the LOW temperature.

How to Commence Melting Later in the Day:

1. To commence melting later in the day, you will need to set your timer.
2. The Timer on your Deluxe Paraffin Spa works like a coffee machine timer. Decide the time that you want to use the Paraffin Spa and then set the hours and minutes to count down to that selected time. For example, if you want to use the Paraffin Spa at 8:00 p.m. and the time is now 12:00 p.m., you will need to set your timer to count down from 8 hours.

How to Set the Timer:

1. Turn the unit on by pressing the power "On/Off" button. When the unit is turned on, the temperature will default to medium, and the word "MED" will show on the display.
2. Press the "TEMP" button once to select the "HIGH" temperature setting. This will ensure that the wax is completely melted when the timer reaches 00:00. **It is recommended that the wax be melted on "HIGH" temperature to achieve the fastest melt time.**
3. Press the "AUTO" key. The display will show 4:00 and flash.
4. Press the "HOUR" key to select the desired number of hours from now that you want to use the Deluxe Paraffin Spa. You can fast forward in 8 digits per second by pressing and holding down the "HOUR" key.
5. Press the "MINUTE" key to select the desired minutes. You can fast forward in 8 digits per second by pressing and holding down the "MINUTE" key.
6. The minimum time setting is 4 hours and the maximum time setting is 23 hours and 59 minutes.
7. When the desired time is set, press the "AUTO" key and the countdown to 00:00 (noted on the display) will commence.
8. When the timer reaches 00:00 on the display, your Deluxe Paraffin Spa is ready to use.
9. Your Deluxe Paraffin Spa comes with a plastic foot rest. Place the foot rest in the unit after

the wax has melted,if desired. **Once wax is melted,do not move unit as spillage may occur!**

10. **NOTE: Paraffin Spa has an auto shut-off feature in timer mode.The unit will automatically shut off within 6 hours of reaching 00:00.After 6 hours,the unit will have to be turned back on to use.**
11. To change the temperature of the Deluxe Paraffin Spa while in Timer countdown mode,press the "AUTO" button and the display will show the temperature level (LOW-MED-HIGH) for 5 seconds. Within the 5 seconds, you can change the temperature of your Deluxe Paraffin Spa by pressing the "TEMP" key until the desired temperature setting shows on the display (LOW, MED, HIGH).The newly selected temperature will continue to flash on the display for a few seconds and then the display will go back to showing the countdown time. **NOTE: TIMER MODE WILL BE ABORTED IF "TEMP"BUTTON IS PRESSED BEFORE "AUTO"BUTTON.**
12. To abort the Timer mode and commence melting immediately, simply press the "TEMP" button once. Then,press the "TEMP" button again until desired melt temperature appears on the display (LOW-MED-HIGH). **NOTE: TIMER MODE WILL BE ABORTED IF "TEMP" BUTTON IS PRESSED BEFORE "AUTO" BUTTON.**

How to Use the Pumice Wand

File clean,dry feet with gentle strokes until rough, hard skin is removed.Use the coarse side for hard, callused skin,the fine side for rough, dry skin.

How to Dip Feet

1. With the unit filled to the "max" line, it takes approximately 8–10 uses (5 dips per foot per use) for the wax to fall below the "min"line. Do not allow wax to fall below the "min" line or over the "max"line.
2. Wash and thoroughly dry feet. **For ultimate softness, use the pumice wand (see "How to use pumice wand" instructions) and/or add your favorite lotion to your feet, hands or elbows and rub in gently prior to dipping. Important:Do not use a heat activated or heat sensitive lotion.**
3. Test the temperature of the wax by dabbing a small amount on the inside of your wrist. If wax seems too hot, simply select a lower temperature by pressing the "TEMP" button until the desired temperature setting (LOW, MED,HIGH) shows on the display. It will take approximately one hour (from HIGH temperature) for Paraffin Spa to reach the MED temperature and 2 hours (from HIGH

temperature) for the Paraffin Spa to reach the LOW temperature.

4. Place the plastic foot rest in the Paraffin Spa if desired.
5. Slowly immerse your foot into the Paraffin Spa chamber. Let the paraffin solidify on your foot. Wait at least three seconds (until the paraffin no longer appears shiny) before redipping.
6. Continue to dip your foot in and out of the wax, following Step 5, to build up layers of paraffin. For the best results, do this at least five times. The more layers, the longer the heat will last.
7. To treat your hands or elbows, follow steps as outlined for feet.For hands, slowly immerse with fingers spread slightly apart.
8. To prolong the heat therapy on your feet, see instructions on "How to use your insulated cotton booties and liners." To prolong the heat therapy on your hands,see instructions on "How to use your insulated cotton mitts and liners (select models).
9. After your treatment is complete,remove the booties or mitts see "How to use your cotton booties or cotton mitts"instructions). **Note: cotton mitts only included in select models.** Remove the plastic liners (see "How to use your cotton booties or cotton mitts" instructions) and discard.Peel the hardened wax from your feet, hands or elbows.Or, use the removal wand included to help remove the wax from your feet. **It is recommended that the hardened wax be discarded due to the buildup of dead skin cells attached to the wax. However, you may put the hardened wax back in the chamber, if desired.**
10. Massage in the oils left on your feet, hands or elbows after the wax is peeled off to get the full invigorating and rejuvenating benefits of your Paraffin Wax Treatment.
11. As soon as you have completed your wax treatment, turn the unit off and unplug it. After unplugging the Paraffin Spa, remove the plastic foot rest. **If the foot rest is not removed, the next time you heat up the Paraffin Spa, the wax will take longer to melt. However, this will not be noticeable in the timer mode as the timer is programmed to compensate for this.** If you discard the used wax after each use, we recommend you clean out the unit completely after 40 uses. If you re-use the wax after every use, we recommend you clean out the unit completely after 20 uses. (See "Cleaning Instructions").
12. If, at this time, cleaning is not necessary, simply allow wax to harden completely and then cover unit and put away until next usage. When you are ready to use the Deluxe Paraffin Spa again, simply remove from storage, plug unit into a polarized outlet of proper voltage, turn unit "on" by following the "How to Commence Melting Immediately" instructions or "How to Commence Melting Later in the Day" and "How to Set Timer"instructions.

HOW TO USE YOUR INSULATED COTTON BOOTIES AND LINERS:

1. Dip your foot in and out of Remington's Paraffin Spa 4-5 times (enough times to form a "soft glove").
2. Put your foot inside a plastic liner (to protect the Insulated Cotton Bootie interior).
3. Then, put your foot inside the Insulated Cotton Bootie to enjoy the soothing effect. The Insulated Cotton Bootie helps retain heat for better results.
4. Repeat steps 1-3 for other foot.

Note: Allow the paraffin to remain in place for 10-15 minutes before removing.

5. After the treatment is complete, remove the booties from your feet. Remove the plastic liners and discard. Peel the hardened wax from your feet, or use the removal wand included to help remove the wax from your feet. **It is recommended that the hardened wax be discarded due to the build-up of dead skin cells attached to the wax. However, you may put the hardened wax back in the Paraffin Spa, if desired.**
6. Massage the oils left on your feet after the wax is peeled off to get the full invigorating and rejuvenating benefits of your Deluxe Paraffin Spa.

How to care for your Insulated Cotton Booties:

- To clean your Insulated Cotton Booties, hand wash with mild soap. Line dry. Do not use washing machine or dryer.

HOW TO USE YOUR INSULATED COTTON MITTS AND LINERS (SELECT MODELS):

1. Dip your hand in and out of Remington's Paraffin Spa 4-5 times (enough times to form a "soft glove").
2. Put your hand inside a plastic liner (to protect the Insulated Cotton Mitt interior).
3. Then, put your hand inside the Insulated Cotton Mitt to enjoy the soothing effect. The Insulated Cotton Mitts help retain heat for better results.
4. Repeat steps 1-3 for other hand.

Note: Allow the paraffin to remain in place for 10-15 minutes before removing.

5. After the treatment is complete, remove the hand mitts from your hands. Remove the plastic liners and discard. Peel the hardened wax from your hands. **It is recommended that the hardened wax be discarded due to the build-up of dead skin cells attached to the wax. However, you may put the hardened wax back in the Paraffin Spa, if desired.**
6. Massage the oils left on your hand after the wax is peeled off to get the full invigorating and rejuvenating benefits of your Deluxe Paraffin Spa.

How to care for your Insulated Cotton Mitts:

- To clean your Insulated Cotton Mitts, hand wash with mild soap. Line dry. Do not use washing machine or dryer.

CLEANING INSTRUCTIONS FOR DELUXE PARAFFIN SPA

1. If you discard the used wax after each use, we recommend you clean out the unit completely every 40 uses. If you re-use the wax after every use, we recommend you clean out the unit completely every 20 uses.
2. Turn unit off.
3. **Before unplugging the Deluxe Paraffin Spa, allow paraffin to remain undisturbed for a few minutes. This will allow any sediment to settle to the chamber bottom.**
4. Unplug unit and let the wax "harden" overnight.
5. Directly before cleaning, plug in the unit and press the power "On/Off" button. Observe the condition of the wax around the entire periphery of the chamber and when the wax is visibly melted around the entire edge (approximately 3-5 minutes) turn the unit off and unplug from electrical outlet.
6. Spread an ample amount of paper towels on a hard, flat, heat-resistant work surface and slowly turn the unit upside down to allow the hardened wax to drop from the chamber onto the paper towels. Or, turn unit upside down over a trash can. To clean the empty chamber, simply wipe with paper towels to remove any sediment.
7. Do not use sharp objects to remove solidified paraffin.
8. Do not use abrasives when cleaning the inside of the Paraffin wax chamber as it will damage the finish.
9. Do not scratch or damage Paraffin wax chamber.
10. To clean the outside of the unit, remove any solidified pieces of paraffin and use a non-abrasive all purpose cleaner.

FULL TWO YEAR WARRANTY

Remington Products Company, L.L.C. warrants this product against any defects that are due to faulty material or workmanship for a two-year period from the original date of consumer purchase. This warranty does not include damage to the product resulting from accident or misuse.

If the product should become defective within the warranty period, we will repair it or elect to replace it free of charge. We will return your product, transportation charges prepaid, provided it is delivered prepaid to: Remington Products Company, L.L.C., 60 Main Street, Bridgeport, CT 06604, Attention: Customer Service Department. In Canada: Remington Products (Canada) Inc., 455 Cochrane Drive, Unit #24, Markham, Ontario L3R 9R4, Attention: Customer Service. Include a letter stating the reason for the return.

REMINGTON SHALL NOT BE RESPONSIBLE FOR ANY INCIDENTAL, SPECIAL OR CONSEQUENTIAL DAMAGES.

This warranty gives you specific legal rights and you may also have other rights which vary from state to state. Some states do not allow the exclusion or limitation of incidental, special or consequential damages.

REMINGTON MONEY BACK OFFER

If within 90 days after you have purchased your Remington Paraffin Spa you are not satisfied and would like a refund, take back the item and sales slip indicating purchase price and date of purchase to the retailer from whom the item was purchased. Remington will reimburse all retailers who accept the item within 90 days from the date of purchase. If you have any questions concerning the money back guarantee, please call 1-800-736-4648 in the U.S. only or 1-905-470-9400 in Canada.

This product will bear one or more compliance marks, Listed by Underwriters Laboratories, Inc. 120 VAC, 60 Hz

Questions or Comments: Please call 1-800-736-4648 in U.S. only or 1-905-470-9400 in Canada.

Visit us at <http://www.remington-products.com>

Remington and the Spa Therapy Collection design are trademarks of Remington Corporation, L.L.C.
© 2001 R.P.C., L.L.C. Rev. 4/01 Job#W01-024 Part# 13929

Patent Pending

REMINGTON®

Better products. Better value. Guaranteed.

Remington Products Company, L.L.C.
60 Main Street Bridgeport, CT 06604

Remington Products (Canada) Inc.
455 Cochrane Drive, Unit #24, Markham, Ontario L3R 9R4.

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Ce produit porte un ou plusieurs sigles et est enregistré par Underwriters Laboratories, Inc. 120 V(CA), 60 Hz