


iWand

User Manual

(English)

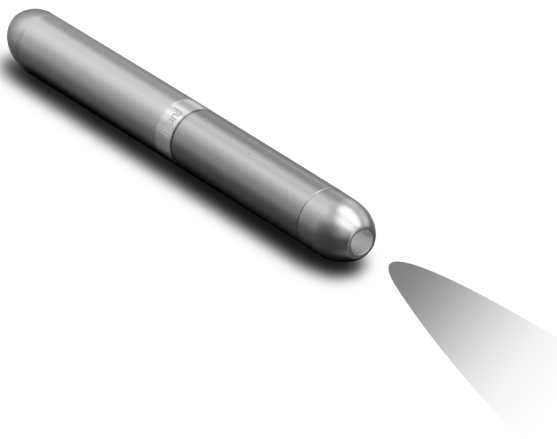

IPC

TABLE OF CONTENTS

Introduction	1
Operation Guide	4
Product Information	7
Warranty	8
Packaging	10
Technical Specifications	10
Disclaimers	11

Cautions

- Do not direct **iWand** toward the eye as it might cause temporary blurred vision.

INTRODUCTION

With the adoption of our PVA field technology, (namely “Proton Vibration Alignment”), **iWand** works by using the medium of light to enhance positive energy flow (known as “Qi” in Chinese medicine) in the body through the skin and acupoint.

Having the advantages of a contactless and non-invasive approach, **iWand** improves microcirculation to achieve the following:

- 1. Rapid absorption of external fluid into skin and deeper tissues*
- 2. Relieve many types of pains and chronic conditions*
- 3. Enhance immunity and promote metabolism for illness prevention*

The benefits of enhancing microcirculation

The standard explanation in the Chinese ancient text is: **Wherever there is blockage of Qi, there is pain.** When Qi flows properly, pain disappears. **iWand** helps to improve Qi flow (or ‘microcirculation’ – see remarks) in the body. The enhanced microcirculation enables rapid absorption of external fluid into the skin and deeper tissues; and improves the physiological functions of the user.

Remarks:

Microcirculation refers to the flow of blood through very thin channels or gaps between cells (smaller than 100 μm). It's responsible for transporting blood cells and enabling the supply of nutrients and chemical substances to/from tissues. In simple words, microcirculation delivers oxygen and nutrients to tissues and carries away metabolic waste from tissues.

Features and Benefits

- Enhance microcirculation



(before use of **iWand**)




(after use of **iWand**)




** Photos show the microcirculation before and after the use of **iWand**. The microcirculation improves after an extended period of use.*

- Rapid absorption of external fluid into skin and deeper tissues
- Relieve many types of pains (headaches, neck and shoulder pain, back pain, knee pain, arthritis, muscle aches, menstrual pain, etc), indigestion, insomnia, skin allergies and other chronic conditions
- Enhance immunity and promote metabolism
- Quick recovery from fatigue and pain caused by poor microcirculation
- Non-invasive Acupuncture
- Applicable for all areas including infected skin, wounds and areas covered with hair
- Free from stimulus response and allergic reaction
- Contactless
- No side effects
- Completely safe and easy to use

OPERATION GUIDE

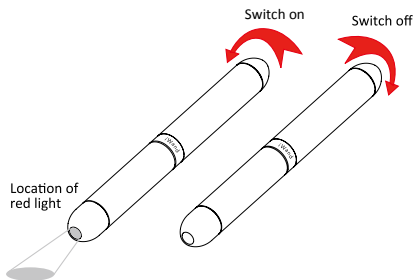


1. Select the targeted skin area. This might be an area covered with external fluid, a pain area or specific acupoint. Turn on the **iWand** and aim the flashing red light onto the selected location. (Suggested usage duration: Absorption purpose - 1 minute; Healthcare purpose - 5 minutes; Pain Relief purpose - 10 minutes)

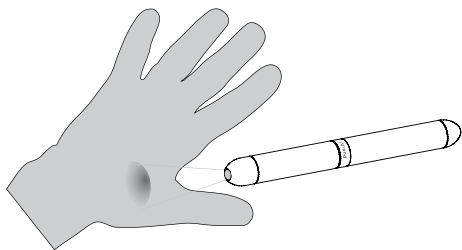

2. When finished, turn off **iWand** until red light is switched off.

Battery Replacement

1. Turn **iWand** in counter clockwise direction to open the casing
2. Insert two 1.5V AAA size batteries (suggested to use Alkaline batteries, do not use rechargeable batteries)
3. Turn **iWand** in clockwise direction to close the casing. Avoid turning **iWand** too tight to turn it on accidentally

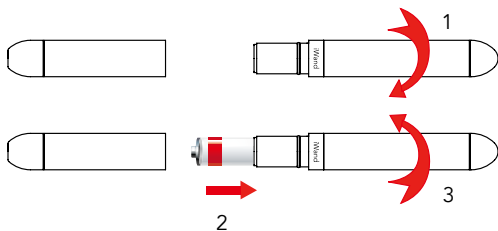

Cautions for Operating iWand

- Do not direct **iWand** toward the eye
- Stop using **iWand** when there is discomfort
- Children or persons who are unable to express themselves, should seek medical advice before use; children must be accompanied by adults when using **iWand**
- For battery operation, **iWand** works only with 1.5V AAA size battery. Do not use rechargeable battery
- Clean **iWand** by dry towel. Do not clean **iWand** with water, detergent, organic solvent or bleach

PRODUCT INFORMATION

Carefully review the following cautions for proper and safe conditions when using **iWand**.

Cautions

- Do not direct **iWand** toward the eye
- Do not use **iWand** in areas of high temperature, high humidity or under direct sunlight
- Do not open, repair, modify, remove parts, make use of components not provided by **iWand**, throw the device from high place or put too much pressure to the device

Storage

- Store **iWand** in the dry environment. Do not place **iWand** in areas of high temperature, high humidity or under direct sunlight
- Remove the batteries if **iWand** will not be used for extended period of time

WARRANTY

There is 1 year warranty service (free labour and components cost) from the purchase date. Customer will be required to present the original invoice stating the serial number for repair service during the warranty period.

Customer will contact the local distributor (details refer to our web site) or dealer for more details about the warranty service. The device must be inspected by **IPC's** appointed authorized representative. The customer will be informed whether the unit can be repaired immediately or longer period is required for repair. **IPC** reserves the right for the final repair period. Local distributor will contact the customer to collect the repaired device after service. If there is no local distributor or repair centre for the customer, he / she should contact the corresponding dealer to arrange the repair service. Customer is responsible for all related courier cost and risk if any.

If there is any additional repair cost involved, the local distributor will inform the customer about the repair cost and repair will take place once confirmed by customer.

If the components or the particular model has been discontinued, the local distributor reserves the right to replace another model with similar function to the customer. The final replaced model and associated cost will be decided by the local distributor.

Importance: This warranty service is excluded from the damage or defective as a result of accident, misuse, abuse, contamination, improper use, unauthorized maintenance, parts removal, damaged or rendered defective caused by incorrect input voltage or using components not provided by the product packaging. If customer chooses the representative not appointed or authorized by **IPC** for repair service, the product warranty will be voided immediately.

PACKAGING

1. **iWand** (1 pc)
2. Portable Pouch (1 pc)
3. User Manual (1 pc)

TECHNICAL SPECIFICATIONS

Size: 15mm(W) x 15mm(H) x 143mm (L)

Weight: 30.4g (battery excluded)

Power: two 1.5V(AAA size) batteries

Maximum power consumption: less than 30 mw

Wavelength of light output: 620 – 645 nanometer

** All provided information is based on the standard testing method designated by Innoworks Product Creation Ltd. They are subject to change without notice.*

DISCLAIMERS

- **IPC** has done the best in preparing this user manual for the overall completeness and correctness. However, **IPC** cannot guarantee any mistake or missing information for the user manual
- **IPC** reserves the right to make changes without any notice on the hardware, software or accessories mentioned in this user manual
- This user manual (or any part of the content) cannot be duplicated, transmitted, copied, stored or translated in any language and/or any format without formal written authorization by **IPC**
- **IPC** assumes no liability, obligation or responsibility for any misuse of **iWand** or other components not provided as part of packaging
- As there are many considerations for the theory of Chinese Acupuncture treatment and meridians, the effectiveness of the **iWand** will be subject to individual application. **IPC** will not guarantee the outcome of individual usage of **iWand**

IPC

Innoworks Product Creation Ltd
www.innoworkspc.com

©2009 Innoworks Product Creation Ltd.