

# SI-200-GIST

## DART AB Sciex Interface Manual User Manual

Manual Revision 2  
Document # 7.5.058



Your source for open-air mass spectrometry

IonSense Inc.  
999 Broadway  
Suite 404  
Saugus, MA 01906

Copyright © 2005-2015 by  
IonSense Inc.  
All rights reserved.

The information in this document has been carefully checked and is believed to be reliable. However, no responsibility is assumed for inaccuracies. Statements in the document not intended to create any warranty, expressed or implied. Specification and performance characteristics of the hardware and software described in the manual may be changed at any time without notice. IonSense Inc. reserves the right to make changes in any product herein in order to improve reliability, design, or function. IonSense does not assume any liability arising out of application or use of any product or circuit described nor does it cover any license under its patent rights or the rights of others.

The apparatus and application of the apparatus described in this document is protected by US Patent Number 6,949,741 and used under license; additional patents pending.

All trademarks are properties of their respective owners.

# Table of Contents

<b>Scope of Delivery</b> .....	4
Standard System Components .....	4
Configuration Dependent Components .....	4
Functional Description .....	5
Safety .....	5
<b>Installation</b> .....	5
Preparing the Source Region .....	5
Attaching the Vapur Interface .....	6
Inserting the Ceramic Ion Transfer Tube into the Vapur Interface .....	7
Attaching the DART SVP source to the Vapur Interface .....	8
Attaching the Mounting Ring (3200, 4000, 5000 only) .....	9
<b>Operation</b> .....	11
Recommended Operating Parameters .....	11
Troubleshooting .....	11
Maintenance .....	11
Replacement Parts .....	12
Revision History .....	12

## Scope of Delivery:

This chapter lists the components of the ABI Sciex Interface.



*Figure 1: Scope of Delivery*

### Standard System Components

- |                       |  |
|-----------------------|--|
| 1. SI-200-A-R         | ABI Sciex Flange Assembly                      |
| 2. SI-999-14          | VAPUR Pump Kit – 1/4"                          |
| 3. JG10019            | Graphite Vespel Ferrule                        |
| 4. SI-200-INT-5500    | Sciex Interlock Plug (installed on SI-200-A-R) |
| 5. SI-200-CC-A        | Sciex Contact Closure Cable                    |
| 6. SHSH_SS188_M6_12mm | M6 x 12 Socket Cap Screws                      |

### Configuration Dependent Components

- |                   |  |
|-------------------|--|
| 1. Varies*        | Ceramic Tubes                            |
| 2. JVLC-ACCUTOF-A | PLATE, SUPPORT, YOKE, ACCUTOF ASSEMBLY** |

**\*NOTE:** See *VAPUR Interface\_Ceramic Tubes Reference Sheet 7.5.038* for ceramic tube lengths (available on [ionsense.com](http://ionsense.com)).

**\*\*NOTE:** Only for 3200, 4000 and 5000 Models.

## Functional Description

The ABI Sciex Interface transfers the ions created by the DART source to the mass spectrometer and allows the DART source to mount on the mass spectrometer.

## Safety

Caution: The surface of the flange can get hot.


## Installation

### Preparing the Source Region

First, before installing the DART interface, place the Sciex Instrument in '**Standby**' mode.

Then, carefully remove all gas and electrical connections from the AB Sciex APCI source and remove the source from the instrument.

The source area should look like the one in the figure below. Carefully remove the cover plate and set that aside. **Caution: the cover plate may be very hot.**

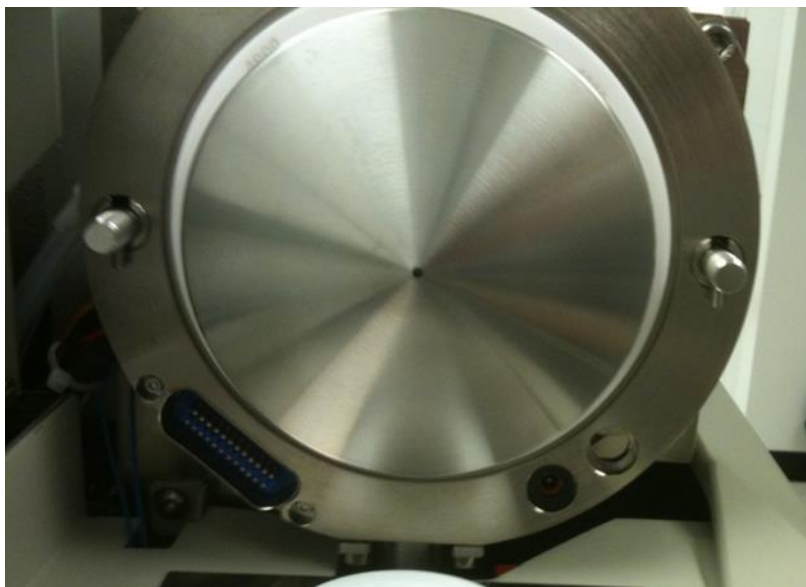

*Figure 2: AB Sciex Source Region*

Once the AB Sciex source and cover plate are removed, the DART flange can be installed.


## Attaching the Vapur Interface

Check that the set screw on the back of the flange is completely inserted into the flange.



**Figure 3:** Set Screw

Orientate the interface so that the turning rods (circled in red in the picture below) line up with the keyholes on either side of the source region. The interlock plug should also line up.

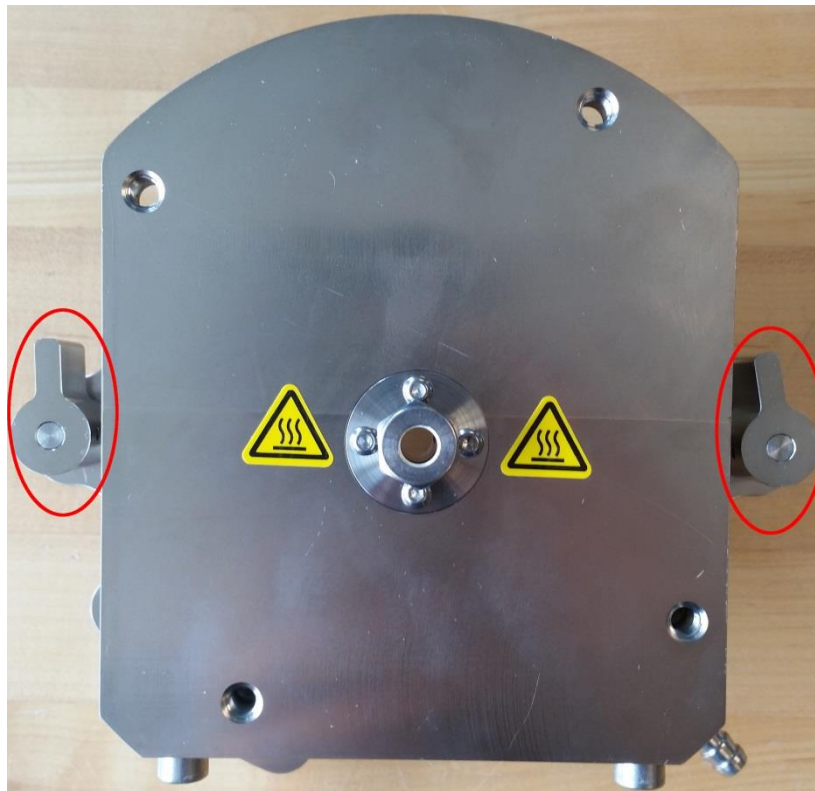

**Figure 4:** Turning Rods


Next gently place the interface over the source region and turn the handles downward locking the interface in place.



**Figure 5:** *Installed Interface*

### **Inserting the Ceramic Ion Transfer Tube into the Vapur Interface**

Remove the Swagelok nut from the front of the Vapur interface and insert the ceramic ion transfer tube. Secure the ceramic tube with one of the provided graphite Vespel ferrules making sure to set a **2 mm** gap between the tube and the APC Inlet.

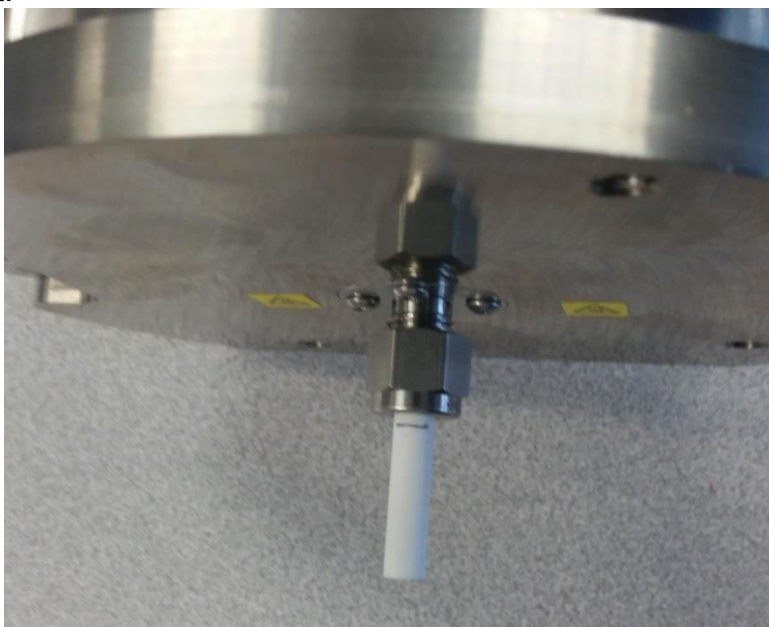

**Figure 6:** *Ceramic ion transfer tube*

**NOTE:** To easily set the 2 mm gap, push the ceramic tube against the Sciex inlet, secure the Swagelok fitting with ferrule and then mark the tube with a pencil. Pull the tube out as shown in Image 3 above so that a 2 mm gap is created inside the interface.


Attach the black rubber pumping line from the diaphragm pump over the barbed pumping port fitting on the side of the Vapur interface.



*Figure 7: Pumping Port*

### **Attaching the DART SVP source to the Vapur Interface**

Attach the base of the source directly to the bottom of the flange with the two M6 x 12 mm socket cap screws provided with the flange.

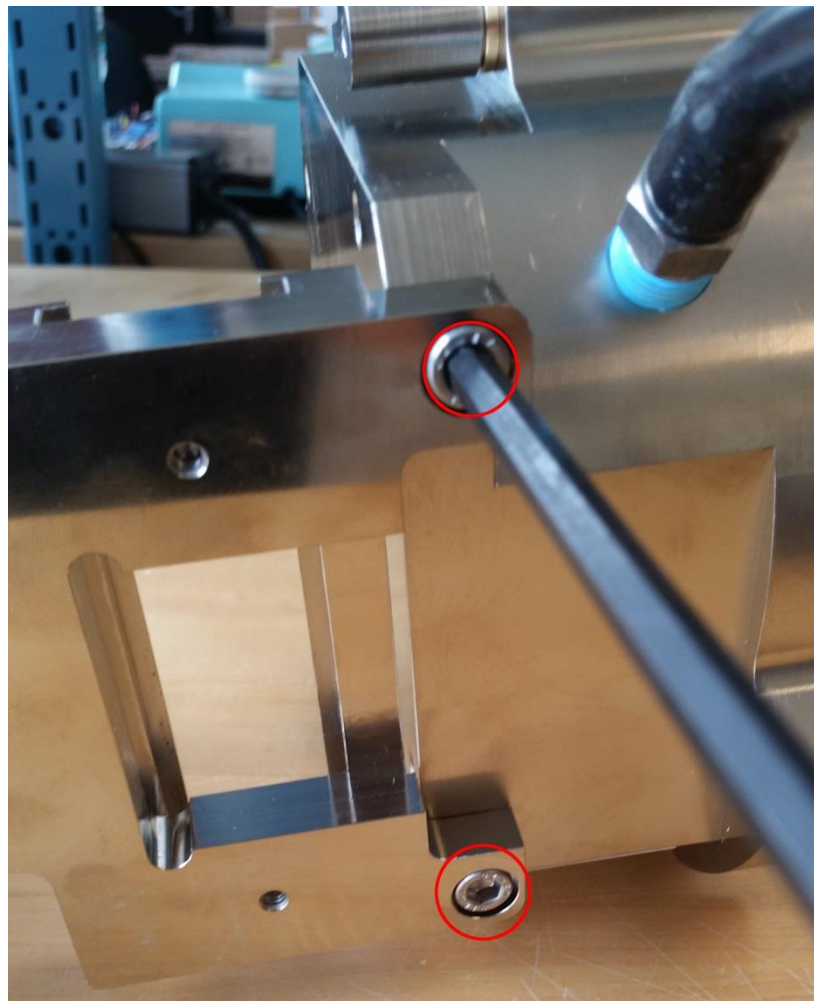


*Figure 8: Attaching the source to the interface*




### Attaching the Mounting Ring (3200, 4000, 5000 only)

Sicex 3200, 4000, 5000 models require a mounting ring before the DART SVP source can be installed.

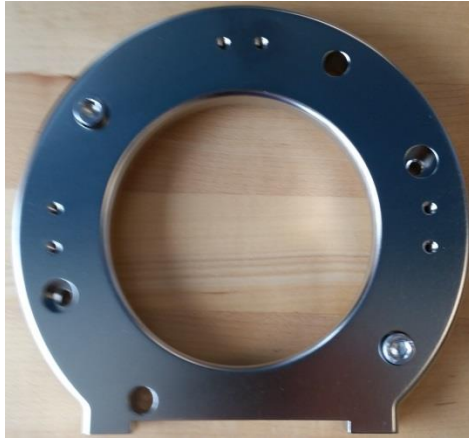


**Figure 9: Mounting Ring**

Install the two alignment screws into the top right and bottom left holes on the flange piece.

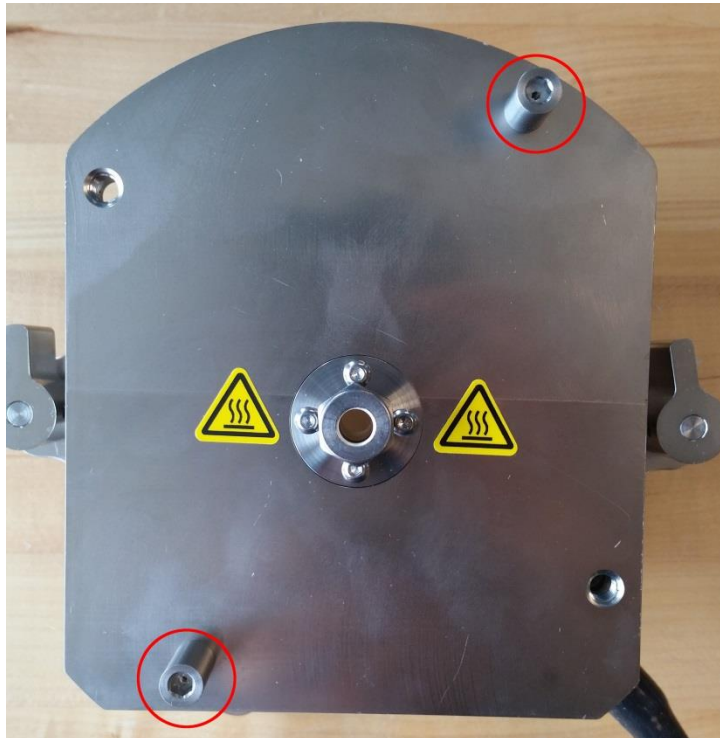


**Figure 10: Alignment Pins**

Use the two M6 x 12 mm socket cap screws to attach the base of the source to the mounting ring.

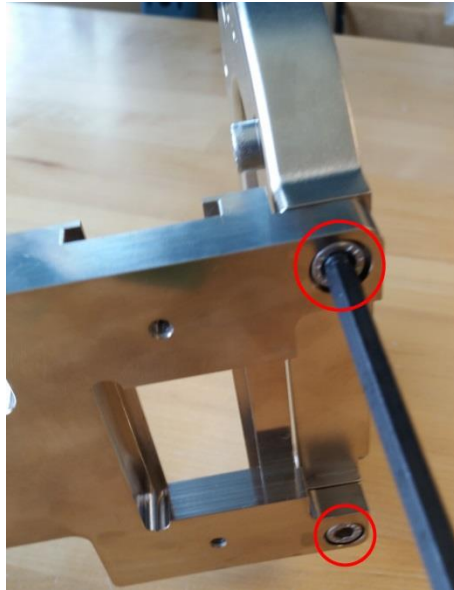

**Figure 11:** Attaching the source to the mounting ring

Slide the mounting ring over the alignment holes and tighten the captive screws on the mounting ring.

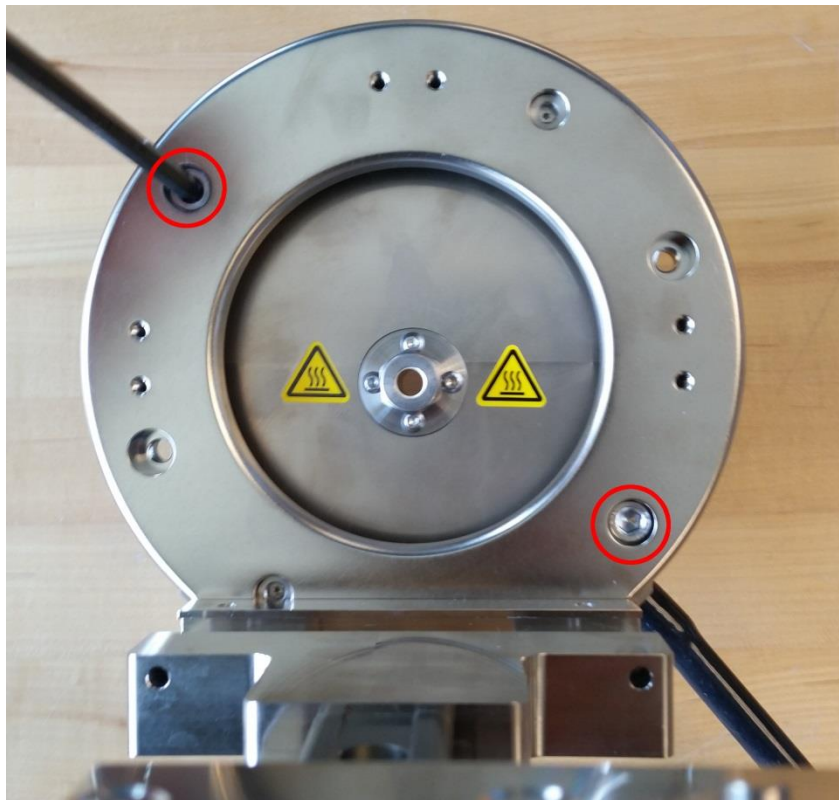

**Figure 12:** Attaching the mounting ring to the interface

# Operation

## Recommended Operating Parameters

Ion Source: Turbo Spray

Curtain Gas (CUR): 10.0

Collision Gas (CAD): ~ 7 (dependent on MRM tune)

IonSpray Voltage (IS) 0.0

Temperature (TEM) 0.0

Ion Source Gas 1 (GS 1) 0.0

Ion Source Gas 2 (GS 2) 0.0

Interface Heater: ON

Suggested Ranges for the “Compounds” Tab (compound-specific):

Declustering Potential (DP): 0-200 V

Entrance Potential (EP): 1-12 V

Collision Energy (CE): 5-130 V

Collision Cell Exit Potential (CXP): 0-58 V

Resolution Q1: Unit

Resolution Q3: Unit

## Troubleshooting

Issue	Solution
Pressure is too high.	<ol style="list-style-type: none"><li>1. Check the vapor pump settings. See the DART Vapor Pump Optimization Protocol 7.5.052 for details (available on ionsense.com).</li><li>2. Check the ceramic tube spacing to ensure a 2 mm gap. See page 8 for further details.</li></ol>

## Maintenance

No special maintenance is required for the ABI Sciex Interface.

## Replacement Parts

SI-999-14	VAPUR Pump Kit - 1/4"
JG10019	Graphite Vespel Ferrule
SI-200-CC-A	Sciex Contact Closure Cable
SI-9xxx-47*	Ceramic Tubes

**\*NOTE:** See VAPUR Interface\_Ceramic Tubes Reference Sheet 7.5.038 for ceramic tube lengths (available on ionsense.com).

## Revision History

Template Revision: 1

Revision History			
REV	DCR #	Description of Change	Effective Date
1	43	Initial Release	11/11/2013
2	221	Updated to new format (User Manual Template_Rev1_7.5.111). Deleted models on page 1. Updated copyright page. Deleted disclaimer sentence above table of contents. Deleted pumping requirements. Deleted and re-wrote pages 5-8 with new pictures.	7/2/2015