

GE Healthcare
Life Sciences

ÄKTA™ avant

Installation Guide

Original instructions

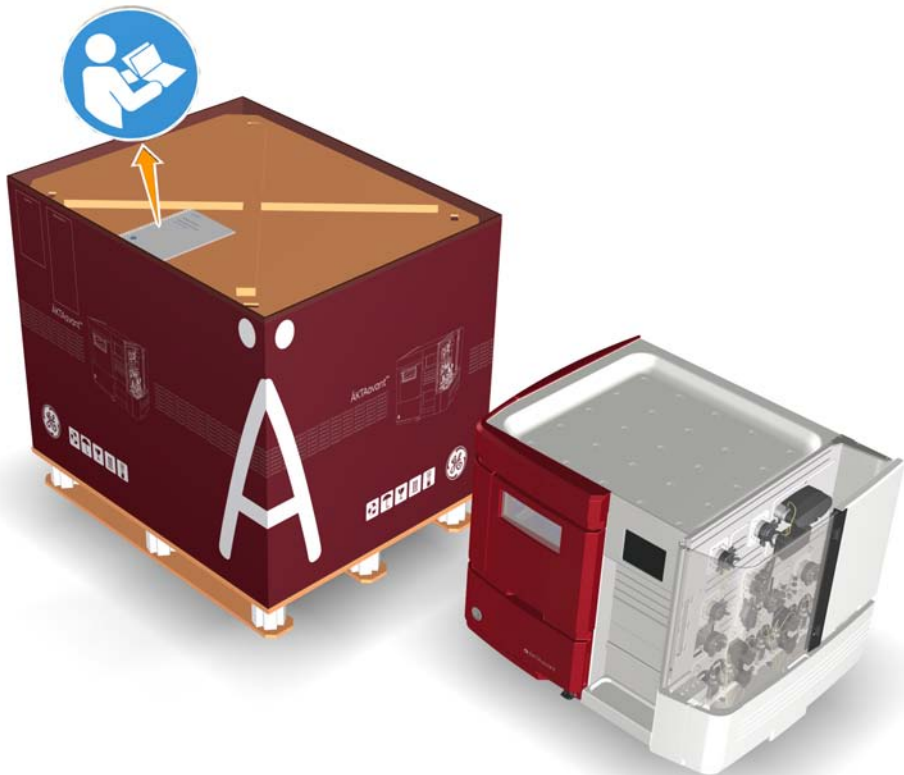


Table of Contents

1	Introduction	5
1.1	Important user information	6
1.2	Regulatory information	8
2	Introduction to ÄKTA avant	12
2.1	The ÄKTA avant instrument	13
2.2	Indicators and controls	14
2.3	Connectors	17
2.4	Network architecture	18
3	Site preparation	20
3.1	Introduction	21
3.2	Delivery and storage	22
3.3	Room requirements	24
3.4	Site environment	28
3.5	Power requirements	29
3.6	Computer requirements	31
3.7	Required materials	32
4	Hardware installation	35
4.1	Unpack the instrument	36
4.2	Install the computer equipment	46
4.3	Connect system units	47
4.4	Install the instrument	50
5	Start the instrument and the computer	57
6	Software installation	58
7	Start UNICORN and connect to system	59
8	Prime inlets and purge pumps	62
8.1	Prime buffer inlets and purge System pump A and System pump B	63
8.2	Prime sample inlets and purge the Sample pump	71
8.3	Prime Q inlets and purge System pump A and System pump B	76
9	Performance tests	81
9.1	Air sensor test	82
9.2	Fraction collector test	85
9.3	Q valve test	90
9.4	System test	96
9.4.1	<i>Prepare the test</i>	97
9.4.2	<i>Run the test</i>	98
9.4.3	<i>Evaluate the test</i>	100

Table of Contents

9.4.4	<i>Possible causes of a failed test</i>	105
9.5	Test protocol for Q valve test	107
9.6	Test protocol for System test	109
10	Reference information	111
10.1	System specifications	112
10.2	Component specifications	114
10.3	Wetted materials	121
10.4	Chemical resistance guide	124

1 Introduction

About the Installation Guide

The purpose of the Installation Guide is to give the necessary instructions to enable users and service personnel to:

- unpack an ÄKTA avant system delivered from the factory
- install the instrument
- install the computer
- install the software
- verify functionality after installation

Read the entire Installation Guide before starting to install the ÄKTA avant system.

In this chapter

This chapter contains the following sections:

Section	See page
1.1 Important user information	6
1.2 Regulatory information	8

1 Introduction

1.1 Important user information

1.1 Important user information

Read this before using the ÄKTA avant system



All users must read the entire *ÄKTA avant Installation Guide and Getting Started with ÄKTA avant*, before installing, operating, or maintaining the instrument.

Do not operate the ÄKTA avant system in any other way than described in the user documentation. Otherwise, you may be exposed to hazards that can lead to personal injury, and you may cause damage to the equipment.

Intended use

ÄKTA avant is a liquid chromatography system intended for process development. The system can be used to screen for optimal choice of columns, media and running parameters to purify selected proteins.

The ÄKTA avant system is intended for research use only, and shall not be used in any clinical procedures, or for diagnostic purposes.

Safety notices

This user documentation contains WARNINGS, CAUTIONS and NOTICES concerning the safe use of the product. See definitions below.

Warnings



WARNING

WARNING indicates a hazardous situation which, if not avoided, could result in death or serious injury. It is important not to proceed until all stated conditions are met and clearly understood.

Cautions



CAUTION

CAUTION indicates a hazardous situation which, if not avoided, could result in minor or moderate injury. It is important not to proceed until all stated conditions are met and clearly understood.

Notices



NOTICE

NOTICE indicates instructions that must be followed to avoid damage to the product or other equipment.

Notes and tips

Note: *A Note is used to indicate information that is important for trouble-free and optimal use of the product.*

Tip: *A tip contains useful information that can improve or optimize your procedures.*

Typographical conventions

Software items are identified in the text by **bold italic** text. A colon separates menu levels, thus **File:Open** refers to the **Open** command in the **File** menu.

Hardware controls, indicators and connections are identified in the text by **bold** text (e.g., **Power** switch).

Text entries that UNICORN™ generates or that the user must type are represented by a monotype typeface (e.g., `\Program Files\GE Healthcare\UNICORN\bin\UNICORN Instrument Server.exe.config`).

1.2 Regulatory information

Introduction

This section describes the directives and standards that are fulfilled by ÄKTA avant.

Manufacturing information

The table below summarizes the required manufacturing information. For further information, see the EC Declaration of Conformity document.

Requirement	Content
Name and address of manufacturer	GE Healthcare Bio-Sciences AB, Björkgatan 30, SE-751 84 Uppsala, Sweden

CE conformity

This product complies with the European directives listed in the table below, by fulfilling the corresponding harmonized standards. For further information, see the EC Declaration of Conformity document.

Directive	Title
2006/42/EC	Machinery Directive (MD)
2006/95/EC	Low Voltage Directive (LVD)
2004/108/EC	ElectroMagnetic Compatibility (EMC) Directive
1999/5/EC	Radio Equipment and Telecommunications Terminal Equipment (R&TTE) Directive

CE marking



The **CE** marking and the corresponding EC Declaration of Conformity, is valid for the instrument when it is:

- used as a stand-alone unit, or

- connected to other GE Healthcare instruments, or
- connected to other products recommended or described in the user documentation, and
- used in the same state as it was delivered from GE Healthcare, except for alterations described in the user documentation.

International standards

This product fulfills the requirements of the following standards:

Standard	Description	Notes
EN ISO 12100	Safety of machinery. General principles for design. Risk assessment and risk reduction.	EN ISO standard is harmonized with EU directive 2006/42/EC
EN 61010-1, IEC 61010-1, UL 61010-1, CAN/CSA C22.2 No. 61010-1	Safety requirements for electrical equipment for measurement, control, and laboratory use.	EN standard is harmonized with EU directive 2006/95/EC
EN 61326-1, IEC 61326-1 (Emission according to CISPR 11, Group 1, class A)	Electrical equipment for measurement, control and laboratory use - EMC requirements	EN standard is harmonized with EU directive 2004/108/EC
EN 301 489-1, EN 301 489-3	Electromagnetic compatibility and Radio spectrum Matters (ERM); ElectroMagnetic Compatibility (EMC) standard for radio equipment and services.	EN standard is harmonized with EU directives 1999/5/EC and 2004/108/EC.
EN 300 330-2	Electromagnetic compatibility and Radio spectrum Matters (ERM); Short Range Devices (SRD); Radio equipment in the frequency range 9 kHz to 25 MHz and inductive loop systems in the frequency range 9 kHz to 30 MHz.	EN standard is harmonized with EU directive 1999/5/EC.

FCC statement

The ÄKTA avant instrument, excluding transmitted radio frequency energy from the Mixer **M9** and UV detector **U9-D**, complies with FCC 47 CFR Part 15b (Federal Communications Commission (FCC) title 47 of the Code of Federal Regulations (CFR), Part 15b, Radio Frequency Devices).

Note: *This equipment has been tested and found to comply with the limits for a Class A digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference when the equipment is operated in a commercial environment. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instruction manual, may cause harmful interference to radio communications. Operation of this equipment in a residential area is likely to cause harmful interference in which case the user will be required to correct the interference at his own expense.*

The modules Mixer **M9** and UV detector **U9-D**, comply with FCC 47 CFR Part 15c.

Note: *This device complies with Part 15c rules. Operation is subject to the following two conditions:*

- *This device may not cause harmful interference, and*
- *This device must accept any interference received, including interference that may cause undesired operation.*

The user is cautioned that any changes or modifications not expressly approved by the manufacturer could void the user's authority to operate the equipment.

Software declaration of conformity

UNICORN 6 is technically compatible with all relevant sections of FDA 21 CFR Part 11.

A part 11-system assessment checklist is available on request through the local GEHC representative.

Regulatory compliance of connected equipment

Any equipment connected to ÄKTA avant should meet the safety requirements of EN 61010-1/IEC 61010-1, or other relevant harmonized standards. Within EU, connected equipment must be CE marked.

Environmental Conformity

Regulation	Title
2011/65/EU	Restriction of Hazardous Substances (RoHS) Directive
2002/96/EC	Waste Electrical and Electronic Equipment (WEEE) Directive
Regulation (EC) No 1907/2006	Registration, Evaluation, Authorization and restriction of Chemicals (REACH)
ACPEIP	Administration on the Control of Pollution Caused by Electronic Information Products, China Restriction of Hazardous Substances (RoHS)

2 Introduction to ÄKTA avant

About this chapter

The ÄKTA avant system comprises the ÄKTA avant instrument, the UNICORN software and accessories. This chapter gives an overview of the ÄKTA avant system and describes the location of indicators and controls on the instrument.

Illustration of the system

The illustration below shows the ÄKTA avant instrument with UNICORN software installed on a computer.



In this chapter

This chapter contains the following sections:

Section	See page
2.1 The ÄKTA avant instrument	13
2.2 Indicators and controls	14
2.3 Connectors	17
2.4 Network architecture	18


2.1 The ÄKTA avant instrument

Introduction

This section provides an overview of the ÄKTA avant instrument.

Illustration of the main parts of the instrument

The illustration below shows the location of the main parts of the instrument.

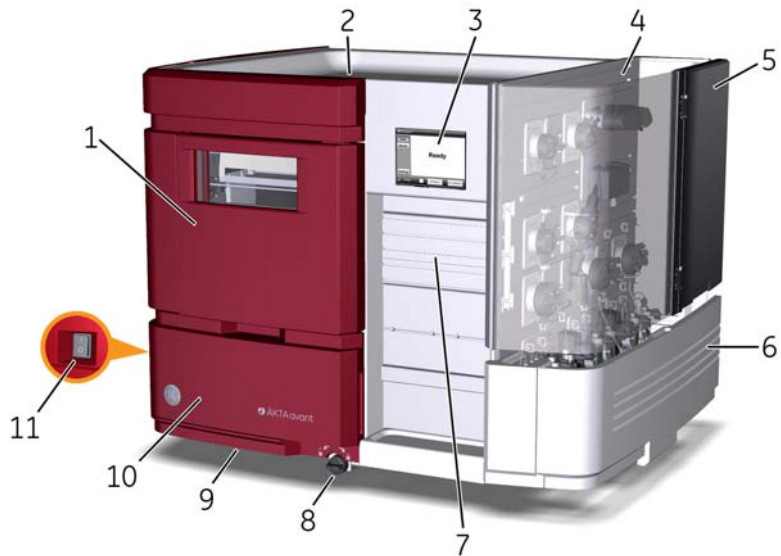

Part	Function	Part	Function
1	Fraction collector	7	Holder rails
2	Buffer tray	8	Swivel foot Lock/Unlock knob
3	Instrument display	9	Swivel foot
4	Wet side	10	Swing out toolbox
5	Foldable door	11	Power switch
6	Pump cover		

2.2 Indicators and controls

Introduction

This section describes the indicators and controls that are available to the user of the ÄKTA avant system.

Illustration

The illustration below shows the location of indicators and controls.



Part	Function
1	Power switch
2	Swivel foot Lock/Unlock knob
3	Instrument display (See description below)




Illustration of the Instrument display

The illustration below shows the Instrument display with the system state **Ready** showing.




Instrument display indicators and buttons

The Instrument display is a touch screen that shows the current system status. The Instrument display includes the following indicators and buttons

Indicator/Button	Description
	Indicates if the Instrument display buttons are unlocked or locked. The buttons can be locked from UNICORN System Control .
	Pauses the run and stops all pumps.
	Resumes instrument operation from the following states: <ul style="list-style-type: none"> • Wash • Pause • Hold

2 Introduction to ÄKTA avant

2.2 Indicators and controls

Indicator/Button	Description
	Indicates that fractionation is ongoing. Do <i>not</i> open the Frac drawer during fractionation.

2.3 Connectors

Introduction

This section describes the connectors for power and communication on the ÄKTA avant instrument.

Illustration

The illustration below shows the location of the connectors.



Part	Function
1	Power input connector
2	Network connector (Ethernet)
3	UniNet-9 connectors Note: <i>Termination plugs must be connected to the connectors that are not in use.</i>

Other connectors are for use by authorized service engineers only.



NOTICE

Misuse of UniNet-9 connectors. The **UniNet-9** connectors at the rear panel should not be mistaken for Firewire connectors. Do not connect any external equipment to the **UniNet-9** connectors. Do not disconnect or move the **UniNet-9** bus cable.

2.4 Network architecture

Introduction


The ÄKTA avant instrument is connected to the UNICORN computer through the Network connector (Ethernet) at the back of the instrument. There are two possible network configurations:

- **Workstation** configuration: local network and a local database
- **Network** configuration: distributed network using TCP/IP, and a central database. Office computers can be used to set up runs, view and evaluate run data.

In this guide, only the workstation configuration is described. For further information about network configuration, please refer to *UNICORN Administration and Technical Manual*.


Illustration, Workstation configuration

The illustration below shows a typical workstation configuration.

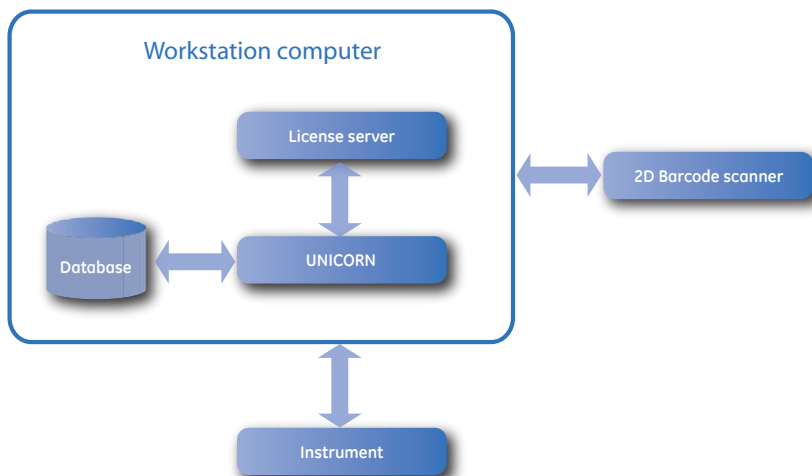

Illustration, Network configuration

The illustration below shows a possible network configuration.

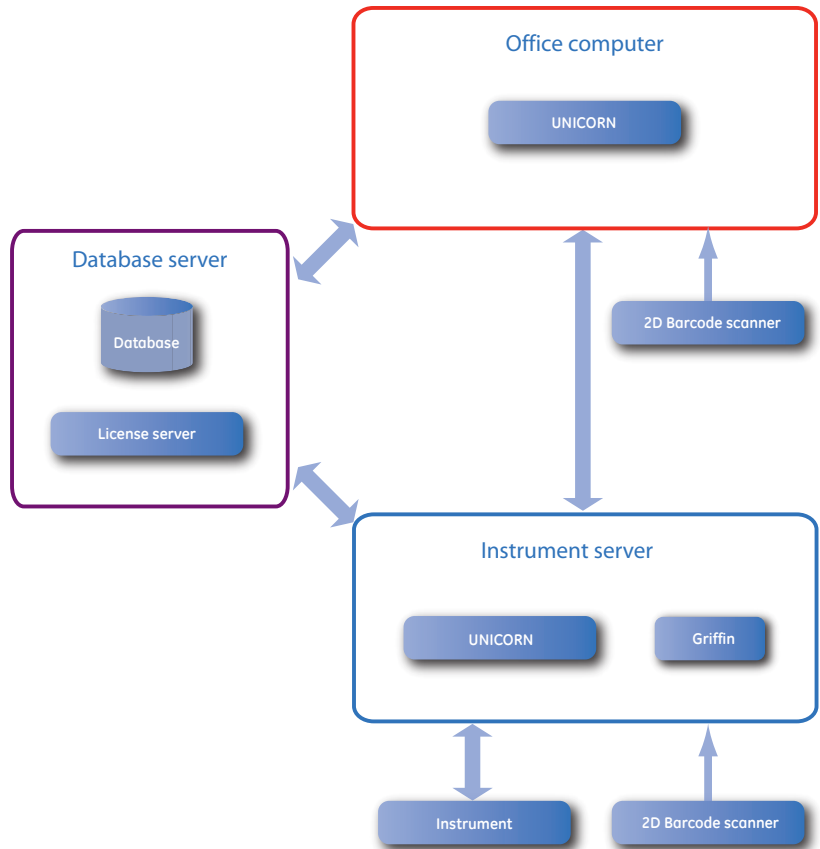

Note: **Griffin** is a software tool that can only be used by GE Healthcare service staff. This service tool is used for diagnostic, testing and quality control.

3 Site preparation

About this chapter

This chapter describes the site planning and the preparations necessary to perform before installation of an ÄKTA avant system. The purpose is to provide planners and technical staff with the data needed to prepare the laboratory for the installation.

In this chapter

This chapter contains the following sections:

Section	See page
3.1 Introduction	21
3.2 Delivery and storage	22
3.3 Room requirements	24
3.4 Site environment	28
3.5 Power requirements	29
3.6 Computer requirements	31
3.7 Required materials	32

3.1 Introduction

Why site preparation?

The laboratory site must be planned and prepared before installing the ÄKTA avant system. The performance specifications of the system can be met only if the laboratory environment fulfills the requirements stated in this chapter. The time spent in preparing the laboratory will contribute to the long term performance of the systems.

Required knowledge

Personnel that will install the system must have:

- A general understanding of how a PC and Microsoft™ Windows™ operating system work. In most cases universal computer functions will not be explained.
 - An understanding of the concepts of liquid chromatography. Terminology and functions will be explained only when they differ from normal practice.
 - A general understanding of risks associated with handling chemicals and substances for use in liquid chromatography.
-

3.2 Delivery and storage

Introduction

This section describes the requirements for receiving the delivery box and storing the instrument before installation.



WARNING

Heavy object. The ÄKTA avant instrument weighs about 116 kg. Use proper lifting equipment, or use four or more persons when moving the instrument. All lifting and moving must be performed in accordance with local regulations.

When you receive the delivery

- Record on the receiving documents if there is any apparent damage on the delivery box. Inform your GE Healthcare representative of such damage.
- Move the delivery box to a protected location indoors.

Delivery box

ÄKTA avant instruments are shipped in a delivery box with the following dimensions and weight:

Contents	Dimensions (mm)	Weight
ÄKTA avant instrument with accessories	w1000 × h900 × d800	155 kg


Storage requirements

The delivery boxes should be stored at a protected place indoors. The following storage requirements must be fulfilled for the unopened boxes:

Parameter	Allowed range
Ambient temperature, storage	-25°C to +60°C
Relative humidity	20% to 95%, non condensing

Equipment for transportation

The following equipment is recommended for handling the delivery boxes:

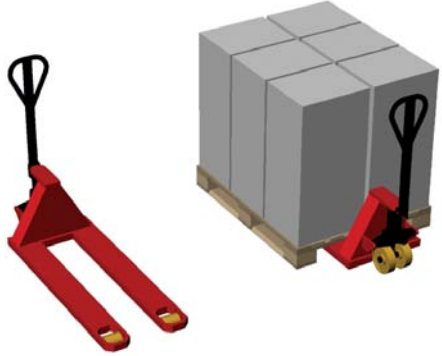Equipment	Specifications
Pallet mover	Suitable for a lightweight pallet 80 × 100 cm 
Cart for transporting the instrument to the lab	Dimensioned to accommodate the size and weight of the instrument

3.3 Room requirements

Introduction

This section describes the requirements for the transportation route and the room where the ÄKTA avant instrument is placed.



WARNING

Protective ground. The ÄKTA avant instrument must always be connected to a grounded power outlet.



WARNING

Only use grounded power cords delivered or approved by GE Healthcare.



WARNING

Do not block access to the power switch and power cord. The power switch must always be easy to access. The power cord with plug must always be easy to disconnect.



WARNING

Explosion hazard. To avoid building up an explosive atmosphere when using flammable liquids, make sure that the room ventilation meets the local requirements.

Transportation route

Doors, corridors and elevators must have a minimum width of 65 cm to allow for transporting the instrument. Allow additional space for moving around corners.

Space requirements

The illustration below shows the space recommended for the ÄKTA avant system.



Allow space on the laboratory bench for:

- handling of samples and buffers (2 × 30 cm)
 - computer and monitor (80 cm)
 - access for service (see below)
-

Service access

To access the rear panel, the instrument can be rotated on a swivel foot. There must be at least 20 cm additional space on the bench to allow for free rotation.



WARNING

Rotating the instrument. Make sure that there is always at least 20 cm of free space around the ÄKTA avant instrument to allow for sufficient ventilation and rotation on the swivel foot. When rotating the instrument, take care not to stretch or squeeze tubing or cables. A disconnected cable may cause power interruption or network interruption. Stretched tubing may cause bottles to fall, resulting in liquid spillage and shattered glass. Squeezed tubing may cause increase in pressure, or block liquid flow. To avoid the risk of knocking over bottles, always place bottles on the buffer tray, and close the doors before rotating the instrument.

Laboratory bench

The bench must be clean, flat and stable to support the weight of the ÄKTA avant system, see table *Equipment weight* below.

Equipment dimensions

The outer dimensions of the ÄKTA avant instrument are shown in the illustration below.



Equipment weight

Item	Weight
ÄKTA avant instrument	116 kg
Computer	9 kg
Monitor	3 kg
<i>Total</i>	<i>128 kg</i>

3.4 Site environment

Introduction

This section describes the environmental requirements for installation of ÄKTA avant.

Room climate

The following requirements must be fulfilled:

- The room must have exhaust ventilation.
- The instrument should not be exposed to direct sunlight.
- Dust in the atmosphere should be kept to a minimum.

Allowed temperature and humidity ranges are specified in the table below.

Parameter	Allowed range
Ambient temperature, operating	4°C to 35°C
Ambient temperature, storage	-25°C to +60°C
Relative humidity, operating	20% to 95%, non condensing

Heat output

The heat output data is listed in the table below.

Component	Heat output
ÄKTA avant instrument	800 W
Computer, incl. monitor and printer	300 W
<i>Total</i>	<i>1100 W</i>

3.5 Power requirements

Introduction

This section describes the power supply requirements for ÄKTA avant.



WARNING

Protective ground. The ÄKTA avant instrument must always be connected to a grounded power outlet.



WARNING

Only use grounded power cords delivered or approved by GE Healthcare.



WARNING

Do not block access to the power switch and power cord. The power switch must always be easy to access. The power cord with plug must always be easy to disconnect.



WARNING

Supply voltage. Make sure that the supply voltage at the wall outlet corresponds to the marking on the instrument, before connecting the power cord.

Requirements

The table below specifies the power requirements.

Parameter	Requirement
Supply voltage	100-240 V ~
Frequency	50-60 Hz
Max power consumption	800 VA

3 Site preparation

3.5 Power requirements

Parameter	Requirement
Number of sockets	1 socket per instrument, up to 3 sockets for computer equipment
Type of sockets	EU or US plugs. Grounded mains sockets, fused or protected by equivalent circuit breaker.
Location of sockets	Maximum 2 m from the instrument (due to length of mains cable). Extension cables can be used if required.

Quality of power

The mains power supply must be stable and conform to specifications at all times to ensure reliable operation of ÄKTA avant. There should be no transient or slow changes in average voltage outside the limits specified above.

3.6 Computer requirements

Introduction

ÄKTA avant systems are controlled by UNICORN software running on a PC. The PC can be part of the delivery or be supplied locally.

The PC used must fulfill the recommendations stated in this section.

Computer specifications

Refer to *UNICORN Administration and Technical Manual* for computer specifications.

3.7 Required materials

Introduction

This section describes the accessories required for the installation and operation of ÄKTA avant.

Buffers and solutions

The buffers and solutions listed below are required during the installation procedure and should be provided at the installation site.

Buffer/solution	Required volume	Scope of use
Distilled water	1 liter	Air sensor test, Fraction collector test, Q valve test, System test
1% acetone in distilled water	0.5 liter	Q valve test
1% acetone and 1 M NaCl in distilled water	0.5 liter	System test
20% ethanol	200 ml	Priming of the pump piston rinsing system

Laboratory equipment

The equipment listed below is required during the installation procedure and should be provided at the installation site.

Equipment	Specification
Flasks, liquid containers	For buffers and waste
Gloves	For protection
Protective glasses	For protection

Fraction collector tubes

The tubes used in the Fraction collector of ÄKTA avant must fulfill the requirements listed in the table below. Examples of manufacturers are also listed in the table.

Tube size [ml]	Diameter [mm]		Height [mm]		Max. volume [ml]	Examples of manufacturers
	Min.	Max.	Min.	Max.		
3	10.5	11.5	50	56	3	NUNC
8	12	13.3	96	102	8	BD Biosciences, VWR
15	16	17	114	120	15	BD Biosciences
50	28	30	110	116	50	BD Biosciences

Deep well plates

Requirements

The deep well plates used in the Fraction collector of ÄKTA avant must fulfill the requirements listed in the table below.

Property	Specification
No. of wells	24, 48, or 96
Shape of wells	Square, not cylindrical
Well volume	10, 5, or 2 ml

Approved deep well plates

The plates listed in the table below are tested and approved by GE Healthcare to be used with ÄKTA avant.

3 Site preparation

3.7 Required materials

Plate type	Manufacturer	Part no.
96 deep well plate	Eppendorf™	951033405/0030 501.306
	BD Biosciences	353966
	Greiner Bio-One	780270
	Porvair Sciences	219009
	Seahorse Bioscience™	S30009
	Whatman™	7701-5200
48 deep well plate	Seahorse Bioscience	S30004
	Whatman	7701-5500
24 deep well plate	Seahorse Bioscience	S30024
	Whatman	7701-5102

4 Hardware installation

About this chapter

This chapter describes the installation procedure of an ÄKTA avant system.

In this chapter

This chapter contains the following sections:

Section	See page
4.1 Unpack the instrument	36
4.2 Install the computer equipment	46
4.3 Connect system units	47
4.4 Install the instrument	50

4.1 Unpack the instrument

Introduction

This section describes how to unpack the ÄKTA avant instrument, and how to lift the instrument onto the bench.



WARNING


Heavy object. The ÄKTA avant instrument weighs about 116 kg. Use proper lifting equipment, or use four or more persons when moving the instrument. All lifting and moving must be performed in accordance with local regulations.

Lift the instrument onto the bench using an overhead crane

It is recommended to use an overhead crane when lifting the ÄKTA avant instrument. Follow the instruction below to lift the instrument onto the bench using an overhead crane and to remove the transport fixations. If an overhead crane is not available, see *Lift the instrument onto the bench by hand*, on page 38.

Step	Action
------	--------

- | | |
|---|------------------------------|
| 1 | Lift off the cardboard hood. |
|---|------------------------------|

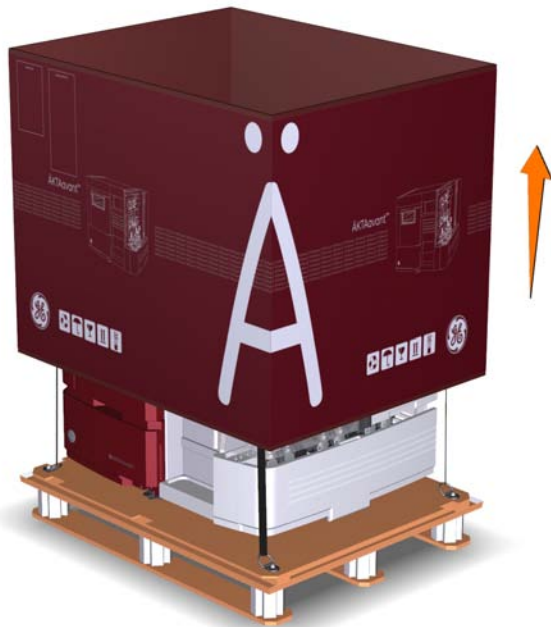


- | | |
|---|--|
| 2 | Make sure that the safety hooks of the straps are properly connected to the lifting yokes located in each of the four corners of the plywood board under the instrument. |
| 3 | Connect both straps to the hook of the overhead crane. |

4 Hardware installation

4.1 Unpack the instrument

Step	Action
------	--------

- | | |
|---|--|
| 4 | Lift the instrument onto the bench using the overhead crane. |
|---|--|

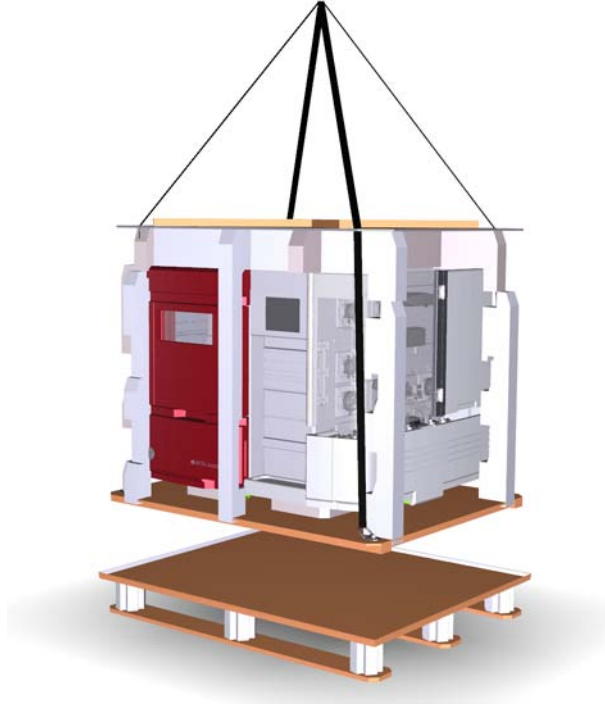

- | | |
|---|--|
| 5 | Disconnect the safety hooks of the straps from the lifting yokes of the plywood board. |
|---|--|


Lift the instrument onto the bench by hand

If an overhead crane is not available, lift off the cardboard hood and disconnect the safety hooks of the straps from the lifting yokes of the plywood board under the instrument. Lift the instrument by hand using the handles illustrated in steps 6-7 in instruction *Remove transport fixations*, on page 39.

Remove transport fixations

Follow the instruction below to remove the transport fixations from the instrument.

Step	Action
1	Lift off the plywood board from the top of the instrument.
2	Remove the pieces of foam from around the instrument.
3	Check the contents in the Buffer tray according to the illustration <i>Accessories packages, on page 42</i> , and lift off the packages from the tray.
4	Unscrew the screws from the Buffer tray securing plate, and remove the plate. Fasten the screws back into place.

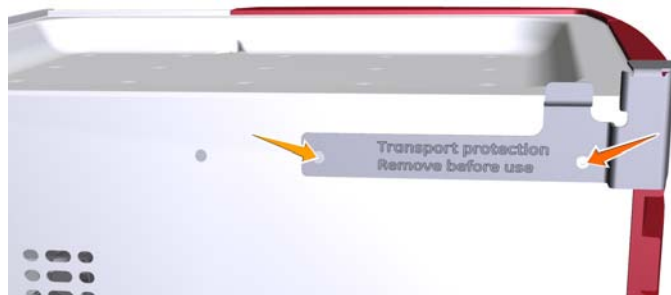


4 Hardware installation


4.1 Unpack the instrument

Step	Action
------	--------



- | | |
|---|---|
| 5 | Remove the pieces of adhesive tape placed on the locations marked in the illustrations below. |
|---|---|

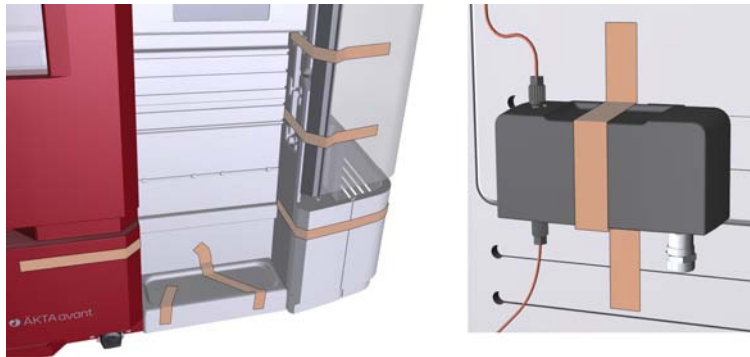

- | | |
|---|--|
| 6 | Open the pump cover, and lift off the tray on the wet side and the tray on the front of the instrument to access the instrument handles. |
|---|--|

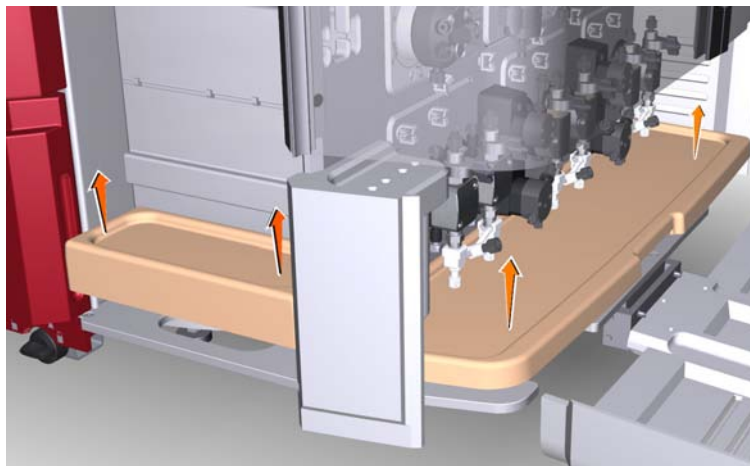


Step	Action
7	Grip the handles on the front and wet side of the instrument. On the back and left-hand side, grip the metal plates located at the bottom of the instrument.

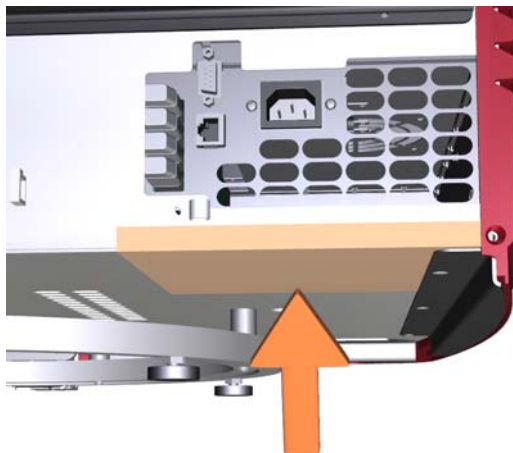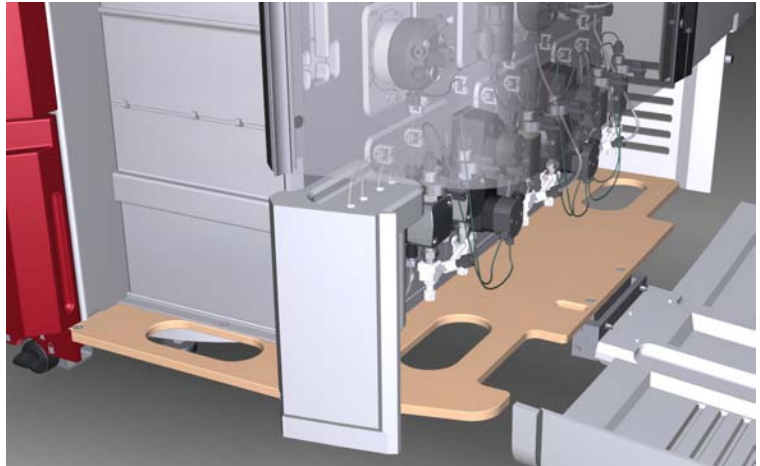

4 Hardware installation
4.1 Unpack the instrument


Step **Action**

- 8 Lift the instrument, and at the same time pull away the plywood board from under the instrument.

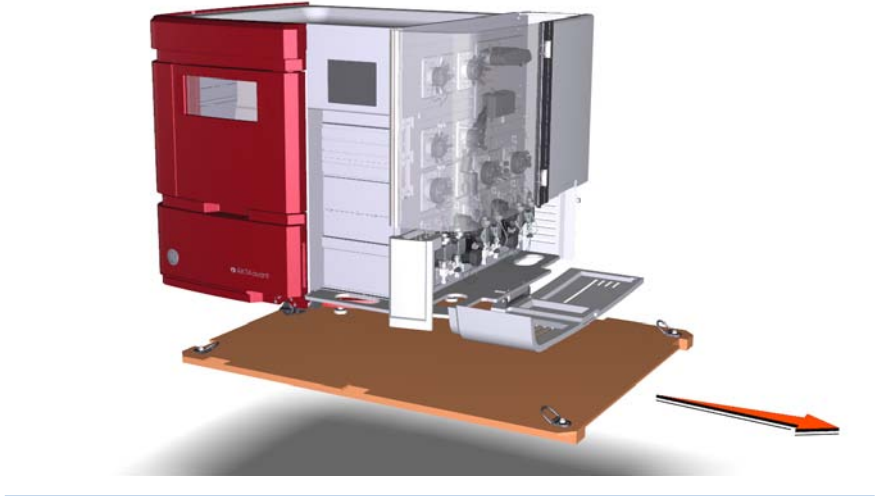

Accessories packages

The illustration below shows the accessories packages placed in the Buffer tray at delivery.

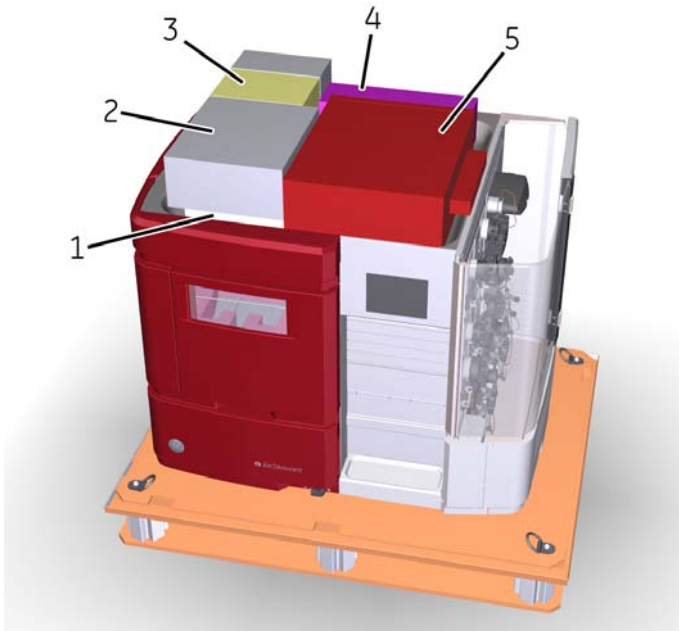


Part	Description
1	Cassette tray
2	Manual box
3	Fraction collector cassettes
4	pH electrode
5	Accessories box

Remove transport fixations from inside the Fraction collector

Follow the instruction below to remove the transport fixations from inside the Fraction collector.

Step	Action
------	--------

- | | |
|---|---|
| 1 | Pull out the deep well plates, the Waste funnel and the Button cover from the holes in the piece of foam located in the Frac chamber. |
|---|---|

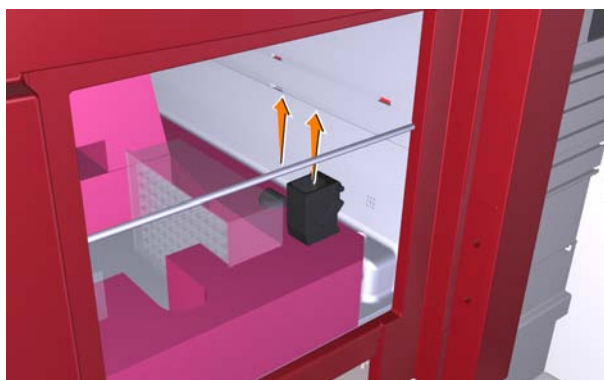


- | | |
|---|--|
| 2 | Fold forward the vertical section of the foam, and fold over the right edge of the foam. |
|---|--|

4 Hardware installation
4.1 Unpack the instrument

Step	Action
------	--------

- | | |
|---|---|
| 3 | Pull out the foam from inside the Fraction collector. |
|---|---|

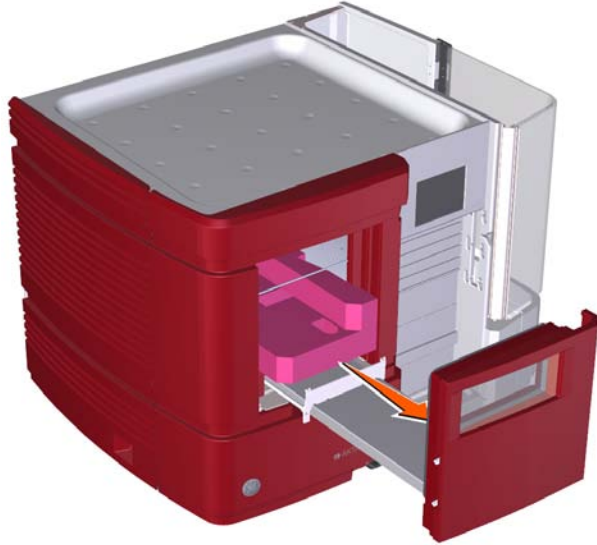


- | | |
|---|--|
| 4 | Pull out the F-shaped piece of foam that holds the Frac arm in position during transportation. |
| 5 | Attach the Waste funnel in position inside the Frac chamber. |



Step	Action
------	--------

- | | |
|---|---|
| 6 | Press the Button cover onto the button on the Frac arm. |
|---|---|

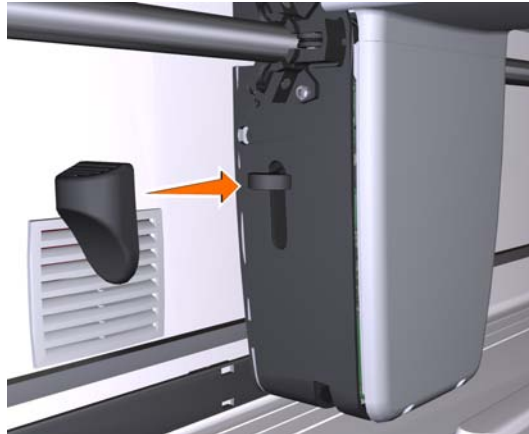

4.2 Install the computer equipment

Introduction

The computer is either:

- Supplied as a part of the ÄKTA avant delivery
 - Supplied locally
-

Unpacking and installing

Unpack and install the computer according to the manufacturer's instructions.

4.3 Connect system units

Introduction

The following interconnections must be made:

- Power supply to the ÄKTA avant instrument
- Power supply to the computer equipment
- Network connection between the computer and the ÄKTA avant instrument.



WARNING

Only use grounded power cords delivered or approved by GE Healthcare.



WARNING

Supply voltage. Make sure that the supply voltage at the wall outlet corresponds to the marking on the instrument, before connecting the power cord.

Connect power to the ÄKTA avant instrument

Follow the instruction below to connect power to the ÄKTA avant instrument.

Step	Action
1	Select the correct power cord to be used. Each instrument is delivered with 2 alternative power cords: <ul style="list-style-type: none">• Power cord with US-plug, 2 m• Power cord with EU-plug, 2 m Discard the unused power cord.
2	Connect the power cord to the Power input connector on the back of the instrument and to a grounded wall outlet 100-240 V ~, 50/60 Hz.

4 Hardware installation

4.3 Connect system units

Step	Action
------	--------

- | | |
|---|---|
| 3 | Attach the power cord to the rear of the instrument using the cable clip. |
|---|---|



Connect power to computer equipment

Follow the manufacturer's instructions to connect power to the:

- computer
- monitor
- local printer, if used

Connect to network

Follow the instructions below to make network connections.

Step	Action
------	--------

- | | |
|---|---|
| 1 | Connect a network cable between the network connector (Ethernet) on the back of the instrument and the computer network card marked for ÅKTA.
The illustration below shows the symbol of the Ethernet connector. |
|---|---|



Step	Action
2	<p data-bbox="435 269 1206 364">If the computer is to be connected to an external network, connect a network cable between the main network card of the computer and a network wall outlet.</p> <p data-bbox="435 382 1206 513">Note: <i>If the computer has not been supplied by GE Healthcare and if network configuration is to be used, see Administration and Technical Manual for further information on network settings.</i></p>

4.4 Install the instrument

Introduction

This section describes how to install the Barcode scanner and the pH electrode, and how to prime the pump piston rinsing systems.

Install the 2D Barcode scanner

Connect the cable of the 2D Barcode scanner to the scanner and to a USB port on the computer.

Install the pH electrode

If pH monitoring is to be used, you need to replace the dummy electrode mounted at delivery with a pH electrode.



CAUTION


pH-electrode. Handle the pH-electrode with care. The glass tip may break and cause injury.

Follow the instruction below to install the pH electrode.

Step	Action
1	Unpack the pH electrode. Make sure that the electrode is not broken or dry.
2	Unscrew the dummy electrode from the flow cell.
3	Pull off the plug from the connector on the front of the pH valve, and store the plug together with the dummy electrode.
4	Remove the cover from the tip of the pH electrode.
5	Carefully insert the electrode in the flow cell. Tighten the locking ring by hand to secure the electrode.
6	Connect the pH electrode cable to the connector on front of the pH valve.

Illustration of the pump piston rinsing systems

The illustration below shows the tubing configuration of the pump piston rinsing systems.

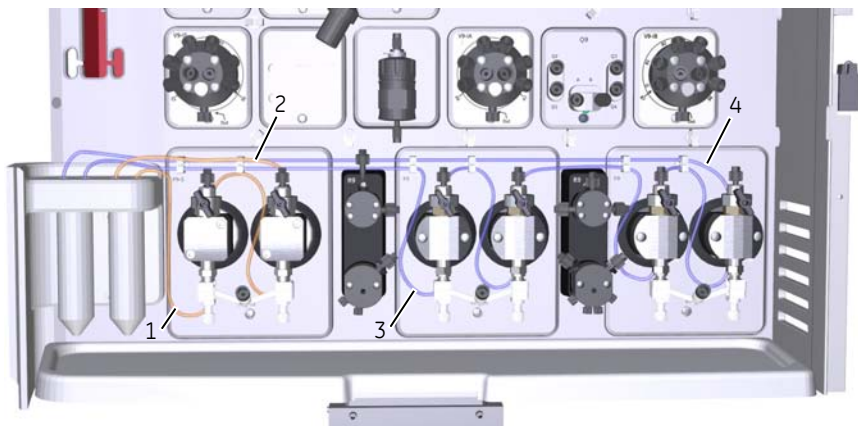

Part	Description
1	Inlet tubing to the Sample pump piston rinsing system
2	Outlet tubing from the Sample pump piston rinsing system
3	Inlet tubing to the System pump piston rinsing system
4	Outlet tubing from the System pump piston rinsing system

Prime the pump piston rinsing system

Follow the instruction below to fill the pump piston rinsing systems with rinsing solution. See the tubing configuration of the rinsing systems in *Illustration of the pump piston rinsing systems*, on page 51.

Step	Action
------	--------

- | | |
|---|--|
| 1 | Unscrew the rinsing system tubes from the holders. |
|---|--|



- | | |
|---|---|
| 2 | Fill each of the rinsing system tubes with 50 ml of 20% ethanol. |
| 3 | Screw the rinsing solution tubes back in the holders. |
| 4 | Immerse the inlet tubing to the System pump piston rinsing system in one of the rinsing solution tubes. |

Note:


Make sure that the inlet tubing reaches to the very bottom of the rinsing solution tube.

- | | |
|---|---|
| 5 | Immerse the inlet tubing to the Sample pump piston rinsing system in the other rinsing solution tube. |
|---|---|

Note:

Make sure that the inlet tubing reaches to the very bottom of the rinsing solution tube.

Step	Action
6	Connect a 25-30 ml syringe to the outlet tubing of the System pump piston rinsing system. Draw liquid slowly into the syringe.

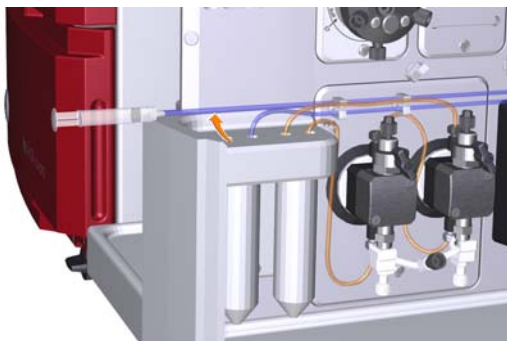

7	Disconnect the syringe and discard its contents.
8	Immerse the outlet tubing in the rinsing solution tube where the inlet tubing of the System pump piston rinsing system is immersed.
9	Connect a 25-30 ml syringe to the outlet tubing from the Sample pump piston rinsing system. Draw liquid slowly into the syringe.
10	Disconnect the syringe and discard its contents.
11	Immerse the outlet tubing in the rinsing solution tube where the inlet tubing of the Sample pump piston rinsing system is immersed.
12	Fill the rinsing solution tubes so that each of the tubes contains 50 ml of 20% ethanol.

Location of waste tubing

All waste tubing is found on the rear of the instrument, see the illustration below.



Part	Description
1	Waste tubing from the Injection valve, the pH valve and the Outlet valve (W , W1 , W2 and W3).
2	Waste tubing from the Fraction collector and the Buffer tray.

Prepare the waste tubing

Follow the instruction below to prepare the waste tubing.

- | Step | Action |
|------|--|
| 1 | Place the four pieces of waste tubing from the Injection valve, the pH valve and the Outlet valve (W , W1 , W2 and W3) in a vessel placed below the bench. |



NOTICE

The maximum level of the waste vessel for the waste tubing from the valves must be lower than 30 cm above the lab bench.

- | | |
|---|--|
| 2 | Place the three pieces of waste tubing from the Fraction collector and the Buffer tray in a waste vessel placed below the bench. |
|---|--|



NOTICE

The maximum level of the waste vessel for the waste tubing from the Fraction collector and the Buffer tray must be lower than the bench height.

- | | |
|---|---|
| 3 | Cut the waste tubing from the Fraction collector and the Buffer tray to appropriate length. It is important that the tubing is not bent and will not be submerged in liquid during the run. |
|---|---|



Note:

If the tubing is too short, replace it with new tubing. Do not lengthen the tubing as this might cause obstruction of the tubing and flooding in the Frac chamber.

4 Hardware installation

4.4 Install the instrument



CAUTION

Make sure that the waste vessels will hold all the produced volume of the run. For ÄKTA avant 25, a suitable waste vessel should typically have a volume of 2 to 10 liters. For ÄKTA avant 150, a waste vessel should have a volume of 40 liters.

5 Start the instrument and the computer

Introduction

This section describes how to start the instrument and the computer.

Instruction

Follow the instructions below to start the instrument and the computer.

- | Step | Action |
|------|--|
| 1 | Switch on the instrument by pressing the Power switch to the I position. |



Result: The instrument starts and the Instrument display states **Not connected**.

- | | |
|---|--|
| 2 | Turn on the computer and monitor according to the manufacturer's instructions. |
|---|--|
-

6 Software installation

Introduction

This section gives an overview of the different UNICORN installation types.

Detailed information about software installation and configuration is available in the *UNICORN Administration and Technical Manual*.

Software installations

UNICORN can be installed as listed below:

- a complete UNICORN installation on a stand-alone workstation (Full installation)
- a UNICORN database and license server (Custom installation)
and
- UNICORN software client and instrument server software on a network client station (Custom installation).

It is also possible to:

- define a system as part of the installation
 - configure E-licenses
 - configure Windows settings necessary for the UNICORN process pictures in a network deployment
 - configure firewall settings, when necessary
 - upgrade UNICORN
 - remove UNICORN installations
and
 - set up a system printer.
-

7 Start UNICORN and connect to system

Introduction

This section describes how to start and log on to UNICORN and how to connect the instrument to UNICORN.

Start UNICORN and log on

Follow the instructions below to start UNICORN and log on to the program.

Step	Action
------	--------

- | | |
|---|---|
| 1 | Double-click the UNICORN icon on the desktop.
<i>Result:</i> The Log On dialog opens. |
|---|---|

Note:

*If there is no connection to the database it is still possible to log on to UNICORN and control a running system. The **Log On** dialog will give the option to start **System Control** without a database. Click **Start System Control** to proceed to the next **Log On** dialog.*

7 Start UNICORN and connect to system

Step	Action
------	--------

2	In the Log On dialog:
---	------------------------------

- select **User Name**.
- and
- enter **Password**.

Note:

*It is also possible to select the **Use Windows Authentication** checkbox and enter a network ID in the **User Name** field.*




- click the **Options** button, and select which UNICORN modules to start.



- click **OK**.

Result: The selected UNICORN modules open.

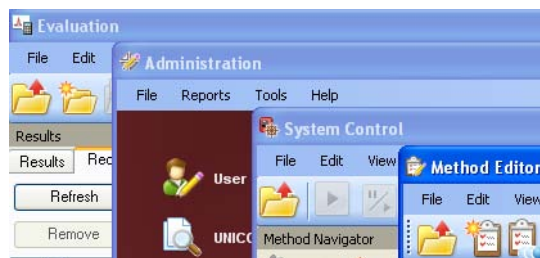

Connect to system


Follow the instructions below to connect the instrument to UNICORN.

Step **Action**

- 1 In the **System Control** module, click the **Connect to Systems** icon.



Result: The **Connect to Systems** dialog opens.

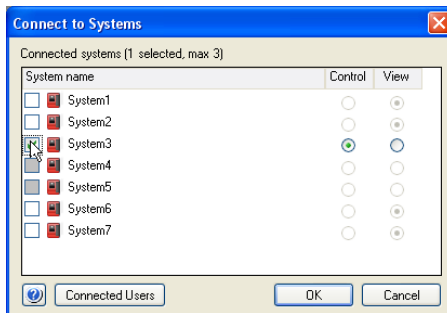

- 2 In the **Connect to Systems** dialog:

- Select a system.
- Select **Control** mode.
- Click **OK**.

Result: The selected instrument can now be controlled by the software.

8 Prime inlets and purge pumps

About this chapter

Before starting the performance tests described in *Chapter 9 Performance tests*, on page 81 it is important to:

- 1 Prime the inlets (fill the inlets with liquid)
- 2 Purge the pumps (remove air from the pumps)

This chapter describes how to prime the buffer inlets, sample inlets, and Q inlets, and how to purge the System pumps and the Sample pump.

In this chapter

This chapter contains the following sections:

Section	See page
8.1 Prime buffer inlets and purge System pump A and System pump B	63
8.2 Prime sample inlets and purge the Sample pump	71
8.3 Prime Q inlets and purge System pump A and System pump B	76

8.1 Prime buffer inlets and purge System pump A and System pump B

Overview

The procedure consists of the following stages:

- 1 Prime all inlet tubing B to be used during the run
 - 2 Prime all inlet tubing A to be used during the run
 - 3 Prepare the system before purging the System pumps
 - 4 Purge System pump B
 - 5 Purge System pump A
 - 6 End the run
-

Prime inlet tubing B

Follow the instructions below to fill all B inlet tubing to be used in the run with appropriate buffer/solution.

Step	Action
1	Make sure that all B inlet tubing that is to be used during the method run is immersed in the correct buffers.


8 Prime inlets and purge pumps

8.1 Prime buffer inlets and purge System pump A and System pump B

Step	Action
------	--------

2 In the **Manual instructions** dialog:


- Select **Flow path:Inlet B**.
- Select the **Position** of the inlet to be filled from the drop-down list. Start at the inlet position with the highest number and end at the inlet position with the lowest number.

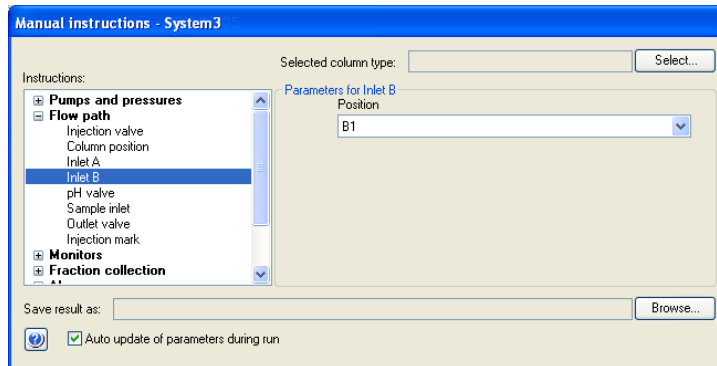

- Click **Execute**.

Result: Inlet valve B switches to the selected port.

3 Connect a 25-30 ml syringe to the purge valve of one of the pump heads of System pump B. Make sure that the syringe fits tightly into the purge connector.

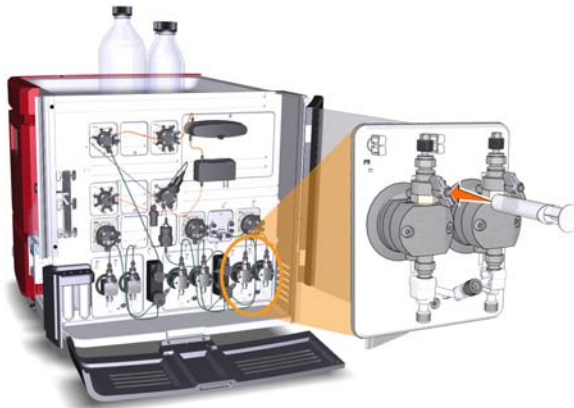


4 Open the purge valve by turning it counter-clockwise about 3 quarters of a turn. Draw liquid slowly into the syringe until the liquid reaches the pump.

Step	Action
5	Close the purge valve by turning it clockwise. Disconnect the syringe and discard its contents.
6	Repeat steps 2-5 for each piece of inlet tubing B that is to be used during the run.

Prime inlet tubing A

Follow the instruction below to fill all A inlet tubing, to be used in the run, with appropriate buffer/solution.


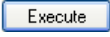
Step	Action
1	Make sure that all A inlet tubing that is to be used during the method run is immersed in the correct buffers.
2	In the Manual instructions dialog: <ul style="list-style-type: none"> • Select Flow path: Inlet A. • Select the Position of the inlet to be filled from the drop-down list. • Click . <p><i>Result:</i> Inlet valve A switches to the selected port.</p>
3	Connect a 25-30 ml syringe to the purge valve of one of the pump heads of System pump A. Make sure that the syringe fits tightly into the purge connector.
4	Open the purge valve by turning it counter-clockwise about 3 quarters of a turn. Draw liquid slowly into the syringe until the liquid reaches the pump.
5	Close the purge valve by turning it clockwise. Disconnect the syringe and discard its contents.
6	Repeat steps 2-5 for each piece of inlet tubing A that is to be used during the run.

8 Prime inlets and purge pumps

8.1 Prime buffer inlets and purge System pump A and System pump B

Prepare the system before purging the System pump heads

Follow the instructions below to prepare the system.

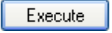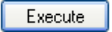
Step	Action
1	Make sure that the piece of waste tubing connected to Injection valve port W1 is placed in a waste vessel.
2	Open the System Control module and select Manual:Execute Manual Instructions . <i>Result:</i> The Manual instructions dialog opens.
3	In the Manual instructions dialog: <ul style="list-style-type: none">• Select Flow path:Injection valve and select System pump waste from the Position drop-down list. Click .<i>Result:</i> The Injection valve switches to waste position. This is necessary to achieve a low backpressure during the purge procedure.• Select Pumps and pressures:System flow. For ÄKTA avant 25, set the Flow rate to 1.0 ml/min. For ÄKTA avant 150, set the Flow rate to 10.0 ml/min. Click .<i>Result:</i> A system flow starts.

Purge System pump B

Follow the instruction below to purge both pump heads of System pump B.


Step	Action
------	--------

- | | |
|---|---|
| 1 | In the Manual instructions dialog: <ul style="list-style-type: none">• Select Pumps and pressures:Gradient.• Set Target to 100% B and Length to 0 min. |
|---|---|

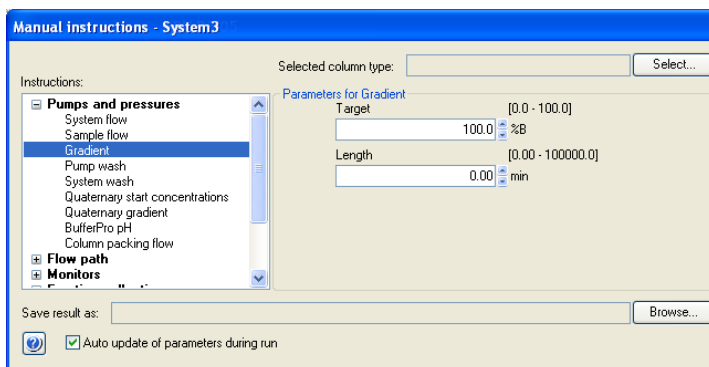

- Click **Execute**.
- Result:* Only System pump B is active.

- | | |
|---|--|
| 2 | In the Manual instructions dialog: <ul style="list-style-type: none">• Select Flow path:Inlet B.• Select the Position of one of the inlets to be used from the drop-down list. |
|---|--|

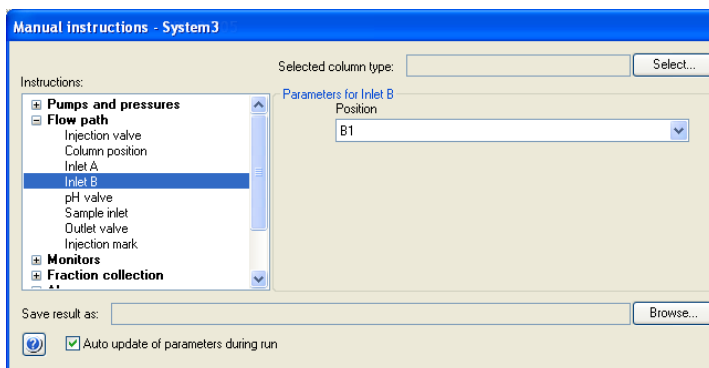


- Click **Execute**.
- Result:* Inlet valve B switches to the selected port.

8 Prime inlets and purge pumps


8.1 Prime buffer inlets and purge System pump A and System pump B

Step	Action
------	--------

- | | |
|---|--|
| 3 | Connect a 25 to 30 ml syringe to the purge valve of the left pump head of System pump B. Make sure that the syringe fits tightly into the purge connector. |
|---|--|

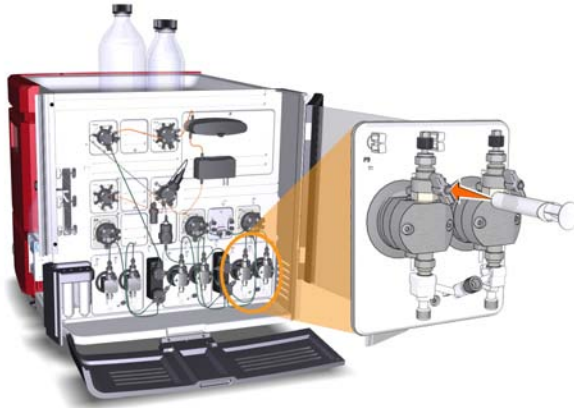


- | | |
|---|--|
| 4 | Open the purge valve by turning it counter-clockwise about 3 quarters of a turn. Draw a small volume of liquid slowly into the syringe (with a rate of about 1 ml per second). |
| 5 | Close the purge valve by turning it clockwise. Disconnect the syringe and discard its contents. |
| 6 | Connect the syringe to the purge valve on the right pump head of System pump B, and repeat steps 4 and 5. Keep the system flow running. |

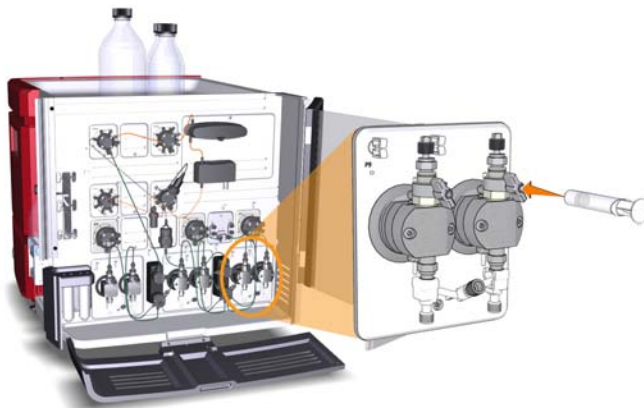

Purge System pump A

Follow the instructions below to purge both pump heads of System pump A.

- | Step | Action |
|------|--|
| 1 | <p>In the Manual instructions dialog:</p> <ul style="list-style-type: none">• Select Pumps and pressures:Gradient.• Set Target to 0% B and Length to 0 min.• Click <input type="button" value="Execute"/>. <p><i>Result:</i> Only System pump A is active.</p> |
| 2 | <p>In the Manual instructions dialog:</p> <ul style="list-style-type: none">• Select Flow path:Inlet A.• Select the Position of one of the inlets to be used from the drop-down list.• Click <input type="button" value="Execute"/>. <p><i>Result:</i> Inlet valve A switches to the selected port.</p> |
| 3 | <p>Connect a 25-30 ml syringe to the purge valve of the left pump head of System pump A. Make sure that the syringe fits tightly into the purge connector.</p> |

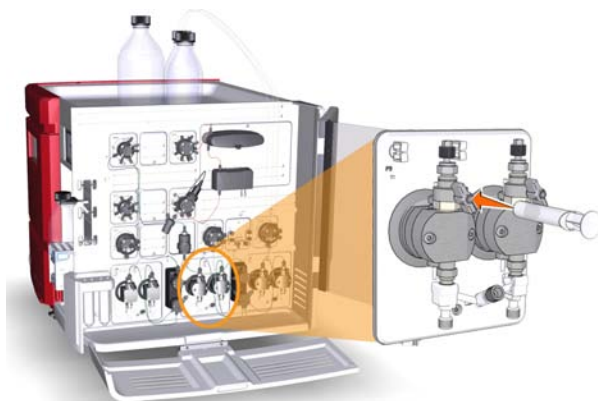


- 4 Open the purge valve by turning it counter-clockwise about 3 quarters of a turn. Draw a small volume of liquid slowly into the syringe (with a rate of about 1 ml per second).

8 Prime inlets and purge pumps

8.1 Prime buffer inlets and purge System pump A and System pump B

Step	Action
5	Close the purge valve by turning it clockwise. Disconnect the syringe and discard its contents.
6	Connect the syringe to the purge valve on the right pump head of System pump A, and repeat steps 5-6. Keep the system flow running.

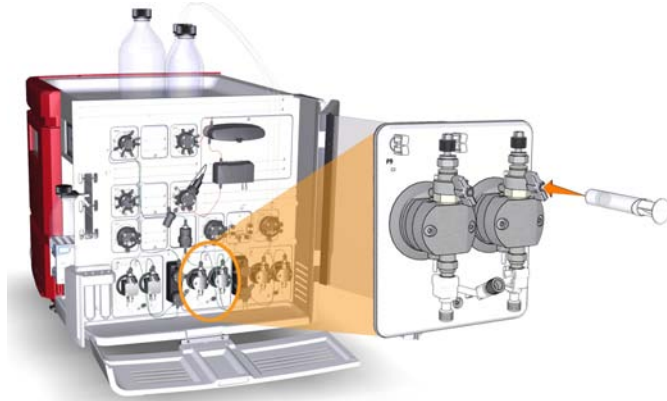

End the run

Click the **End** icon in the **System Control** toolbar to end the run.



8.2 Prime sample inlets and purge the Sample pump

Overview


The procedure consists of the following steps:

- 1 Prime all sample inlet tubing to be used during the run
- 2 Prepare the system before purging the Sample pump
- 3 Purge the Sample pump
- 4 End the run

Prime sample inlet

Follow the instructions below to fill all sample inlet tubing, to be used in the run, with appropriate buffer/solution.

Step	Action
1	Make sure that all sample inlet tubing that is to be used during the method run is immersed in the correct buffers.
2	In the Manual instructions dialog: <ul style="list-style-type: none">• Select Flow path:Sample inlet.• Select the Position of the inlet to be filled from the drop-down list.

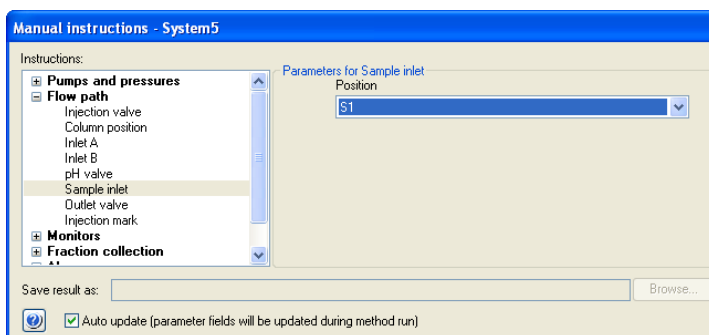

- Click **Execute**.


Result: The Sample valve switches to the selected port.

8 Prime inlets and purge pumps

8.2 Prime sample inlets and purge the Sample pump

Step	Action
------	--------

- | | |
|---|--|
| 3 | Connect a 25-30 ml syringe to one the purge valves of the pump heads of the Sample pump. Make sure that the syringe fits tightly into the purge connector. |
|---|--|

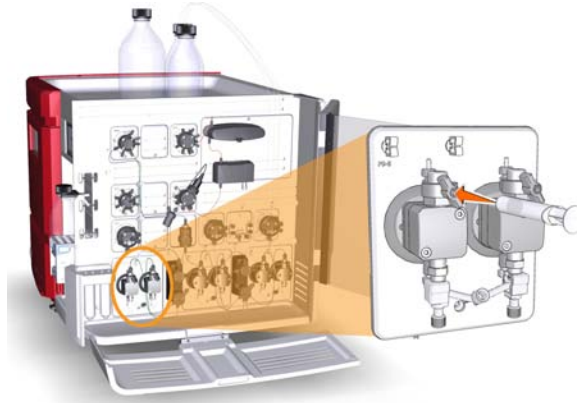



- | | |
|---|---|
| 4 | Open the purge valve by turning it counter-clockwise about 3 quarters of a turn. Draw liquid slowly into the syringe until the liquid reaches the pump. |
| 5 | Close the purge valve by turning it clockwise. Disconnect the syringe and discard its contents. |
| 6 | Repeat steps 2-5 for each sample inlet that is to be used in the method run. |

Prepare the system before purging the Sample pump

Follow the instruction below to prepare the system.


Step	Action
------	--------

- | | |
|---|--|
| 1 | Make sure that the piece of waste tubing connected to Injection valve port W1 is placed in a waste vessel. |
| 2 | Open the System Control module, and select Manual:Execute Manual Instructions .
<i>Result:</i> The Manual instructions dialog opens. |

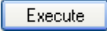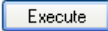Step	Action
3	<p>In the Manual instructions dialog:</p> <ul style="list-style-type: none">• Select Flow path:Injection valve and select Sample pump waste from the Position drop-down list. Click . <p><i>Result:</i> The Injection valve switches to waste position. This is necessary to achieve a low backpressure during the purge procedure.</p> <ul style="list-style-type: none">• Select Pumps and pressures:Sample flow. For ÄKTA avant 25, set the Flow rate to 1.0 ml/min. For ÄKTA avant 150, set the Flow rate to 10.0 ml/min. Click . <p><i>Result:</i> A sample flow starts.</p>

Purge the Sample pump

Follow the instruction below to purge both the pump heads of the Sample pump.

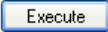
Step	Action
1	<p>In the Manual instructions dialog:</p> <ul style="list-style-type: none">• Select Flow path:Sample inlet.• Select Buffer from the Position drop-down list.• Click . <p><i>Result:</i> The Sample inlet valve switches to Buffer port.</p>
2	<p>Connect a 25-30 ml syringe to the purge valve of the left pump head of the Sample pump. Make sure that the syringe fits tightly into the purge connector.</p>

8 Prime inlets and purge pumps


8.2 Prime sample inlets and purge the Sample pump

Step	Action
------	--------

- | | |
|---|--|
| 3 | Open the purge valve by turning it counter-clockwise about 3 quarters of a turn. Draw a small volume of liquid slowly into the syringe (with a rate of about 1 ml per second). |
|---|--|

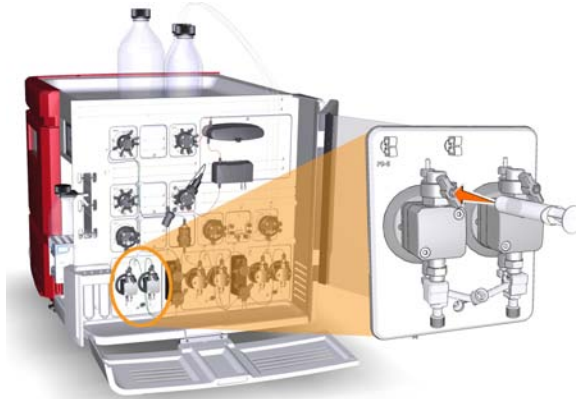

- | | |
|---|---|
| 4 | Close the purge valve by turning it clockwise. Disconnect the syringe and discard its contents. |
| 5 | Connect the syringe to the purge valve on the right pump head of the Sample pump, and repeat steps 2-3. Keep the sample flow running. |

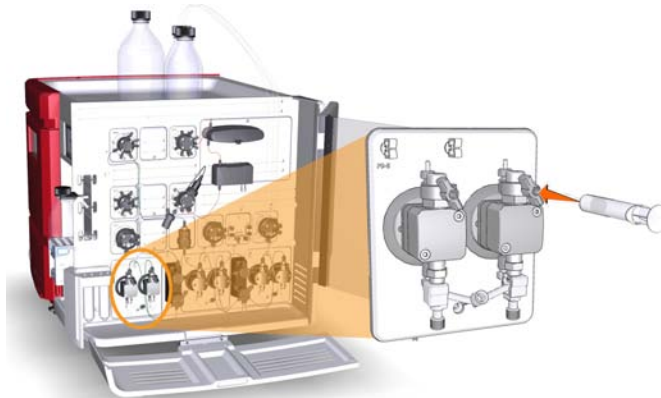

End the run

Click the **End** icon in the **System Control** toolbar to end the run.



8.3 Prime Q inlets and purge System pump A and System pump B

Overview

The procedure consists of the following steps:

- 1 Prepare the system before priming the Q inlets and purging the System pumps
 - 2 Prime all Q inlet tubing
 - 3 Purge the Quaternary valve and the System pumps
 - 4 End the run
-

Prepare the system before priming the Q inlets and purging the System pumps

Note: *Inlet **A1** and **B1** must be immersed in buffer or water. When the System pumps are synchronized, the Inlet valves are positioned to **A1** and **B1** for a short moment.*


Step	Action
------	--------

1	Make sure that the piece of waste tubing connected to Injection valve port W1 of the Injection valve is placed in a waste vessel.
---	--

2	Open the System Control module, and select Manual:Execute Manual Instructions .
---	---

Result: The **Manual instructions** dialog opens.

3	In the Manual instructions dialog:
---	---

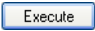- Select **Flow path:Injection valve** and
- choose **System pump waste** from the **Position** drop-down list.
- Click 

Result: The injection valve switches to the waste position.

Note:


This is necessary in order to achieve a low backpressure during the purge procedure.

Step	Action
4	<p>Select Pumps and pressures: System flow:</p> <ul style="list-style-type: none">For ÄKTA avant 25: Set the Flow rate to 1.0 ml/min.For ÄKTA avant 150: Set the Flow rate to 10.0 ml/min. <p>Click <input type="button" value="Execute"/></p> <p><i>Result:</i> A system flow starts.</p>

Prime the Q inlets

Follow the instruction below to prime the Q inlets.

Step	Action
1	Make sure that the pieces of inlet tubing marked A1 , B1 and Q1-Q4 are immersed in the correct buffers.
2	<p>In the Manual instructions dialog:</p> <ul style="list-style-type: none">Select Pumps and pressures: Quaternary start concentrations.Set Start concentration Q1 to 100%. Make sure that the other start concentrations are set to 0%.

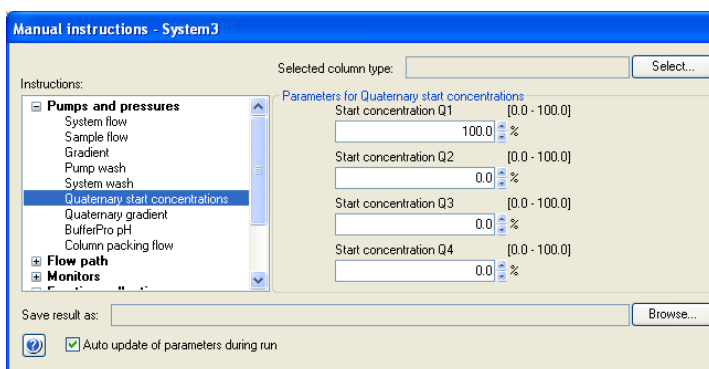


- Click

8 Prime inlets and purge pumps

8.3 Prime Q inlets and purge System pump A and System pump B

Step	Action
------	--------

- | | |
|---|--|
| 3 | Connect a 25-30 ml syringe to one of the purge valves of either of the system pumps. Make sure that the syringe fits tightly into the purge connector. |
|---|--|

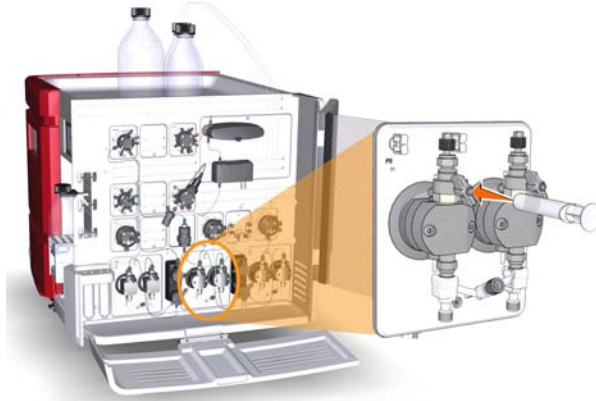

- | | |
|---|--|
| 4 | Open the purge valve by turning it counterclockwise about 3 quarters of a turn. Draw 10 ml of liquid into the syringe. Check that the Q1 inlet is filled with liquid. |
| 5 | Close the purge valve by turning it clockwise. Disconnect the syringe and discard its contents. |
| 6 | Repeat steps 2-5 for Q2 , Q3 and Q4 respectively by setting the respective Quaternary start concentration to 100%. |

Tip:


The piece of inlet tubing that is immersed in distilled water should be the last piece of inlet tubing to be primed.

Purge the Quaternary valve and the System pumps

Follow the instruction below to purge the Quaternary valve, System pump A and System pump B. Note that both pump heads of each System pump have to be purged.

Step	Action
------	--------

- 1 In the **Manual instructions** dialog:
 - Select **Pumps and pressures: Pump wash**, and select **All** from the **BufferPro / Q inlets** drop-down list.

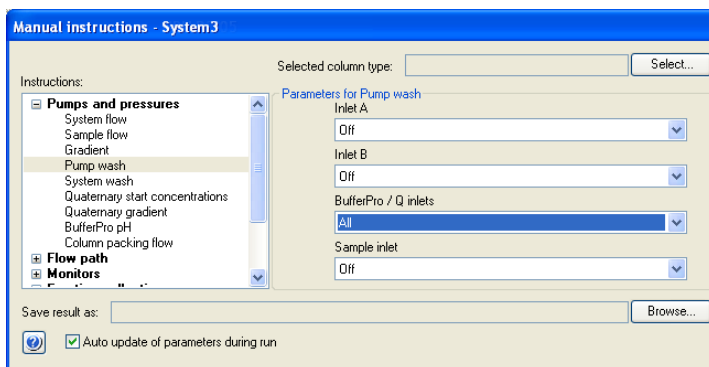

- Click **Execute**.

Result: A simultaneous pump wash of all the Q inlets is started. This will remove air from the Quaternary valve.


- 2 Wait until the pump wash is completed.

8 Prime inlets and purge pumps

8.3 Prime Q inlets and purge System pump A and System pump B

Step	Action
------	--------

- | | |
|---|---|
| 3 | Connect a 25-30 ml syringe to the left purge valve of the selected System pump. Make sure that the syringe fits tightly into the purge connector. |
|---|---|

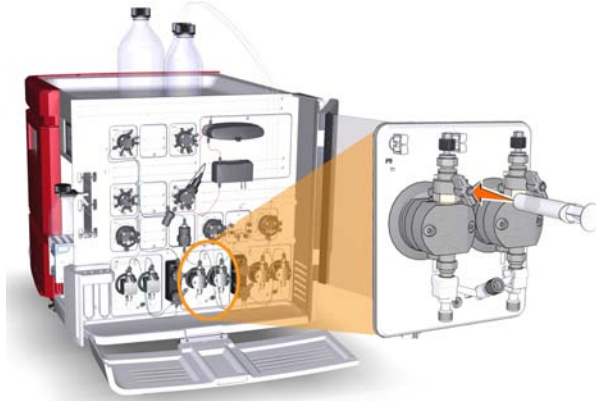

- | | |
|---|--|
| 4 | Open the purge valve by turning it counterclockwise about 3 quarters of a turn. Draw 10 ml of liquid slowly into the syringe with a rate of about 1 ml per second. |
| 5 | Close the purge valve by turning it clockwise. Disconnect the syringe and discard its contents. |
| 6 | Repeat steps 3-5 for the other three purge valves of the System pumps to get rid of air in all pump heads. Keep the system flow running. |

End the run

Click the **End** icon in the **System Control** toolbar to end the run.



9 Performance tests

About this chapter

Performance tests should be run after installation to check the function of the ÄKTA avant system. Performance tests can also be used at any time to check the condition of the system, for example, after a prolonged stop. This chapter describes how to prepare, run, and evaluate the four performance tests available.

In this chapter

This chapter contains the following sections:

Section	See page
9.1 Air sensor test	82
9.2 Fraction collector test	85
9.3 Q valve test	90
9.4 System test	96
9.5 Test protocol for Q valve test	107
9.6 Test protocol for System test	109

9.1 Air sensor test

Method description

The Air sensors A, B, and S are integrated in Inlet valve A, Inlet valve B, and Sample inlet valve. The Air sensor test checks if the air sensors detect air. The method run takes approximately 5 minutes.

Required material

The following material is required:

- Syringe, 25-30 ml
 - Distilled water
-

Prepare the test

Follow the instruction below to prepare the system before method start.

Step	Action
1	Immerse the pieces of inlet tubing marked A1 , B1 , and the piece of sample inlet tubing marked Buff in distilled water.
2	Prime the inlets A1 , B1 , and the sample inlet Buff , and purge the pumps. See <i>Section 8.1 Prime buffer inlets and purge System pump A and System pump B, on page 63</i> and <i>Section 8.2 Prime sample inlets and purge the Sample pump, on page 71</i> .
3	Disconnect the inlet tubing connected to the inlet valve positions A4 , B4 , and S4 . During the test method air is introduced into the inlet valves through these inlet ports to test the function of the air sensors.

Start the test

Follow the instruction below to start the Air sensor test.

Step	Action
1	In the System Control module, select System:Performance Test and Report . <i>Result:</i> The System Performance Test and Report dialog opens.

Step	Action
2	<p>In the System Performance Test and Report dialog:</p> <ul style="list-style-type: none"> • Select Air Sensor Test • Click the Run Performance Method button <p><i>Result:</i> The Start Protocol of the Air sensor test opens.</p>
3	<p>Click Next in the dialogs of the Start Protocol to open the next dialog. The dialogs are described in table <i>Overview of the Start Protocol dialogs, on page 83</i>.</p>
4	<p>In the last dialog of the Start Protocol, Result name and location, click Finish.</p> <p><i>Result:</i> The Air sensor test starts.</p>

Overview of the Start Protocol dialogs

The table below describes the dialogs of the Start Protocol.

Dialog	Description
Notes	Displays the Method Notes of the method. The Method Notes contains a method description and instructions on how to run the method. This dialog also allows the user to enter Start Notes .
Evaluation Procedures	Allows the user to select to save the report to file (recommended) and/or to print the report.
Result Name and Location	Allows the user to change result name and result location.

During the run

A **Message** dialog opens during the run. Read the messages in the dialog, and make sure that necessary preparations have been performed. Click the **Confirm and Continue** button in the **Message** dialog to change system state from **System Pause** to **Run** and proceed with the test. Alternatively, click the **Confirm** button in the **Message** dialog and click the **Continue** button on the Instrument display.

Automatic evaluation

The system automatically generates a report when the test is finished. The report can be printed in two ways:

- If the option **Save the report to file** was selected in the **Evaluation Procedures** dialog in the **Start Protocol** of the test, the report is saved in location `C:\Program Files\GE Healthcare\UNICORN\Temp`. Open the report and print it using your default Windows printer. (Recommended)

Note: *The exact search path depends on the location of the UNICORN installation folder.*

- If the option **Print report** was selected in the **Evaluation Procedures** dialog in the **Start Protocol** of the test, the report is automatically printed on the system printer. Note that a system printer must have been set up, see *UNICORN Administration and Technical Manual*.

Print the report and check the status of the tests. For each of the tests the report states "The test passed" or "The test failed".

Possible causes of a failed test

The table below describes possible causes of a failed test. When possible sources of error have been checked and taken care of, run the test once again.

Cause	Action
Faulty air sensor	For further information, see <i>User Manual, chapter Troubleshooting</i> .
Incorrect preparation of tubing	Make sure that the tubing was correctly prepared, see <i>Prepare the test, on page 82</i> .

9.2 Fraction collector test

Method description

The Fraction collector test checks the functionality of the Fraction collector. The method run takes approximately 10 minutes. For ÄKTA avant 25, fractionation of 2 ml is performed in three sequential wells in each of two 96 deep well plates. For ÄKTA avant 150, fractionation of 20 ml is performed in three sequential tubes placed in each of two Cassettes for 50 ml tubes.

Required material

The following material is required:

- Distilled water
 - Syringe, 25-30 ml
 - For ÄKTA avant 25:
 - Cassettes for deep well plate, 2 pcs
 - 96 deep well plates, 2 pcs
 - For ÄKTA avant 150:
 - Cassettes for 50 ml tubes, 2 pcs
 - 50 ml tubes, 12 pcs
-


Prepare the test

Follow the instruction below to prepare the system before method start.

Step	Action
1	Immerse the piece of inlet tubing marked A1 in distilled water.
2	Prime inlet A1 and purge System pump A. See <i>Section 8.1 Prime buffer inlets and purge System pump A and System pump B, on page 63.</i>
3	Place the deep well plates or tubes in the cassettes. Make sure that the deep well plates are rotated so that the well marked A1 is positioned above the A1 marking on the Cassette.


Step **Action**

- 4 Place the Cassettes on the Cassette tray, one on Cassette position 1 and one on Cassette position 6. Make sure that the Cassette type codes (see illustration below) faces the front of the tray marked with the GE logo.

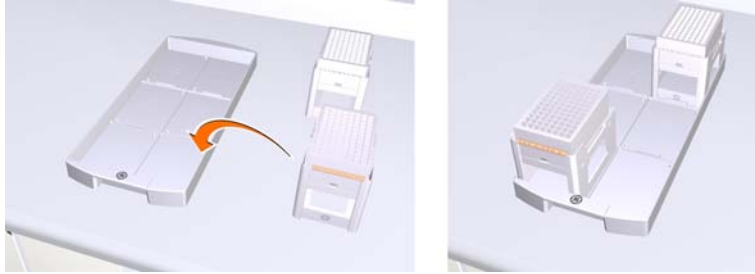

- 5 Open the Frac drawer by pressing the handle upwards, and pulling out the drawer.

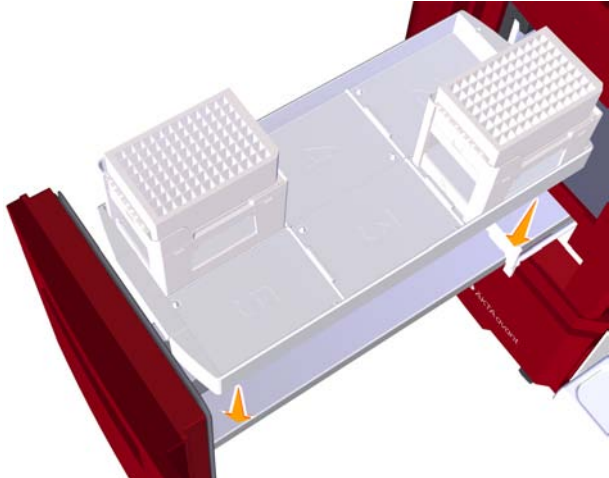

Step	Action
6	Place the Cassette tray on the Tray support of the Frac drawer. Make sure that the front of the tray (marked with the GE logo) faces the front of the drawer.
	
7	Close the Frac drawer. Make sure that it snaps into closed position.
8	In the System Control module, select View:Fraction collection content . In the Fraction collector content dialog, check that the automatic scanning procedures have been performed correctly.
9	In the System Control module, select System:Settings . <i>Result:</i> The System Settings dialog opens.
10	In the System Settings dialog: <ul style="list-style-type: none"><li data-bbox="439 1197 969 1224">• Select Fraction collector:Fractionation settings.<li data-bbox="439 1252 1166 1306">• For ÄKTA avant 25: In the Fractionation mode drop-down list, select Accumulator.<li data-bbox="439 1324 1105 1352">• In the Fractionation order drop-down list, select Row-by-row.<li data-bbox="439 1370 564 1397">• Click OK.

Start the test

Follow the instruction below to start the Fraction collector test.

Step	Action
1	In the System Control module, select System:Performance Test and Report . <i>Result:</i> The System Performance Test and Report dialog opens.
2	In the System Performance Test and Report dialog: <ul style="list-style-type: none">• Select Fraction Collector Test• Click the Run Performance Method button <i>Result:</i> The Start Protocol of the Fraction collector test opens.
3	Click Next in the dialogs of the Start Protocol to open the next dialog. The dialogs are described in table <i>Overview of the Start Protocol dialogs, on page 88</i> .
4	In the last dialog of the Start Protocol, Result name and location , click Finish . <i>Result:</i> The Fraction collector test starts.

Overview of the Start Protocol dialogs

The table below describes the dialogs of the Start Protocol.

Dialog	Description
Notes	Displays the Method Notes of the method. The Method Notes contains a method description and instructions on how to run the method. This dialog also allows the user to enter Start Notes .
Result Name and Location	Allows the user to change result name and result location.

During the run

Check that an accumulator wash is performed at the beginning of the run. During the accumulator wash the Dispenser head is positioned over the Waste funnel and liquid is squirted from the nozzle several times.

Evaluate the result

For ÄKTA avant 25: Check that correct volume, 2 ml, is collected in each of the wells **1A1-1A3** and **6A1-6A3**. Also, check that the fractionation marks in the chromatogram correspond to the filled wells.

For ÄKTA avant 150: Check that correct volume, 20 ml, is collected in each of the tubes **1A1-1A3** and **6A1-6A3**. Also, check that the fractionation marks in the chromatogram correspond to the filled tubes.

Possible causes of a failed test

The table below describes possible causes of a failed test. When possible sources of error have been checked and taken care of, run the test once again.

Cause	Action
<p>Wrong volumes collected in the wells or tubes, and disturbances on the system pressure curves:</p> <ul style="list-style-type: none"> - Air trapped in System pump A - Faulty System pump A 	<p>Air in pumps: Make sure to prime inlet tubing A1 and purge System pump A before method start, see <i>Section 8.1 Prime buffer inlets and purge System pump A and System pump B, on page 63</i>.</p> <p>Faulty pump: See <i>User Manual, chapter Troubleshooting</i>.</p>
<p>Liquid is collected in the wrong wells or tubes:</p> <ul style="list-style-type: none"> - Incorrect fractionation settings 	<p>Check the fractionation settings in the System Settings dialog.</p>
<p>Liquid is collected in wrong Cassette:</p> <ul style="list-style-type: none"> - Incorrect preparation of the Fraction collector 	<p>Make sure that only two Cassettes are placed on the Cassette tray, and that the Cassettes are positioned on the cassette positions 1 and 6.</p>
<p>Fractionation interrupted:</p> <ul style="list-style-type: none"> - Impurities on the Dispenser head of the Fraction collector 	<p>Clean the Fraction collector. See <i>User Manual, chapter Maintenance</i>.</p>
<p>Incorrect preparation of buffer and tubing</p>	<p>Make sure that the system was correctly prepared, see <i>Prepare the test, on page 85</i>.</p>

9.3 Q valve test

Method description

The Q valve test checks the functionality of the Quaternary valve. Correct gradient formation is tested by producing a series of gradient steps. The method run takes approximately 30 minutes for ÄKTA avant 25 and 70 minutes for ÄKTA avant 150.

Required material

The following material is required:

- 1% acetone in distilled water
 - Distilled water
 - Capillary loop **Ref 1**
 - Syringe, 25-30 ml
 - For ÄKTA avant 25: Mixer, 1.4 ml
 - For ÄKTA avant 150: Mixer, 5 ml
-

Prepare the test

Follow the instructions below to prepare the system before method start.

Step	Action
1	Immerse the pieces of inlet tubing marked Q1 and Q3 in 1% acetone in distilled water.
2	Immerse the pieces of inlet tubing marked Q2 and Q4 in distilled water.
3	Connect the capillary loop marked Ref 1 between Column valve ports 1A and 1B to generate a backpressure.
4	For ÄKTA avant 25: Make sure that the Mixer with a chamber volume of 1.4 ml is installed. For ÄKTA avant 150: Make sure that the Mixer with a chamber volume of 5 ml is installed. For further information, see <i>User Manual chapter Maintenance</i> .
5	Prime all Q inlets and purge the System pumps. See <i>Section 8.3 Prime Q inlets and purge System pump A and System pump B, on page 76</i> .

Start the test

Follow the instruction below to start the Q valve test.

Step	Action
1	In the System Control module, select System:Performance Test and Report . <i>Result:</i> The System Performance Test and Report dialog opens.
2	In the System Performance Test and Report dialog: <ul style="list-style-type: none"> • Select Q Valve Test. • Click the Run Performance Method button. <i>Result:</i> The Start Protocol of the Q valve test opens.
3	Click Next in the dialogs of the Start Protocol to open the next dialog. The dialogs are described in table <i>Overview of the Start Protocol, on page 91</i> .
4	In the last dialog of the Start Protocol, Result name and location , click Finish . <i>Result:</i> The Q valve test starts.

Overview of the Start Protocol

The table below describes the dialogs of the Start Protocol.

Dialog	Description
Notes	Displays the Method Notes of the method. The Method Notes contains a method description and instructions on how to run the method. This dialog also allows the user to enter Start Notes .
Evaluation Procedures	Allows the user to select to save the report to file (recommended) and/or to print the report.
Result Name and Location	Allows the user to change result name and result location.

During the run

A **Message** dialog opens during the run. Read the messages in the dialog, and make sure that necessary preparations have been performed. Click the **Confirm and Continue** button in the **Message** dialog to change system state from **System Pause** to **Run** and proceed with the test. Alternatively, click the **Confirm** button in the **Message** dialog and click the **Continue** button on the Instrument display.

Automatic evaluation

The system automatically generates a report when the test is finished. The report can be printed in two ways:

- If the option **Save the report to file** was selected in the **Evaluation Procedures** dialog in the **Start Protocol** of the test, the report is saved in location `C:\Program Files\GE Healthcare\UNICORN\Temp`. Open the report and print it using your default Windows printer. (Recommended)

Note: *The exact search path depends on the location of the UNICORN installation folder.*

- If the option **Print report** was selected in the **Evaluation Procedures** dialog in the **Start Protocol** of the test, the report is automatically printed on the system printer. Note that a system printer must have been set up, see *UNICORN Administration and Technical Manual*.

Print the report and check the status of the tests. For each of the tests the report states "The test passed" or "The test failed".

Manual evaluation - Preparation

If you prefer to evaluate your results manually, follow the instructions below.

Step	Action
------	--------

- | | |
|---|--|
| 1 | In the Evaluation module in UNICORN, click the Open Result Navigator icon. |
|---|--|



Result: The **Result Navigator** opens.

- | | |
|---|---|
| 2 | In the Result Navigator , double-click the test result name to open the test result. |
|---|---|


Step **Action**

- 3 In the **Chromatogram** pane of the **Evaluation** module:
- Select the curve UV1_265.
 - Right-click and select **Vertical marker**.
- 4 Click the vertical marker and drag it to the constant section of each plateau of the curve UV1_265. Read the corresponding absorbance values in the upper left corner of the **Chromatogram** pane. Enter the values in column 2 in the tables in *Section 9.5 Test protocol for Q valve test, on page 107* as follows.

Enter the absorbance value at...	in table...
0% UV (100% Q2)	<i>Step response test result - Q1, Q2</i> <i>Step response test result - Q3, Q4</i> <i>Step response test result - Q1, Q2, Q3, Q4</i>
100% UV (100% Q1)	<i>Step response test result - Q1, Q2</i> <i>Step response test result - Q3, Q4</i> <i>Step response test result - Q1, Q2, Q3, Q4</i>
5% Q1, 95% Q2	<i>Step response test result - Q1, Q2</i>
95% Q1, 5% Q2	<i>Step response test result - Q1, Q2</i>
5% Q3, 95% Q4	<i>Step response test result - Q3, Q4</i>
95% Q3, 5% Q4	<i>Step response test result - Q3, Q4</i>
50% UV (25% Q1, 25% Q2, 25% Q3, 25% Q4)	<i>Step response test result - Q1, Q2, Q3, Q4</i>

Illustration of chromatogram

The illustration below shows a chromatogram from a Q valve test. The gradient steps are marked in the illustration.

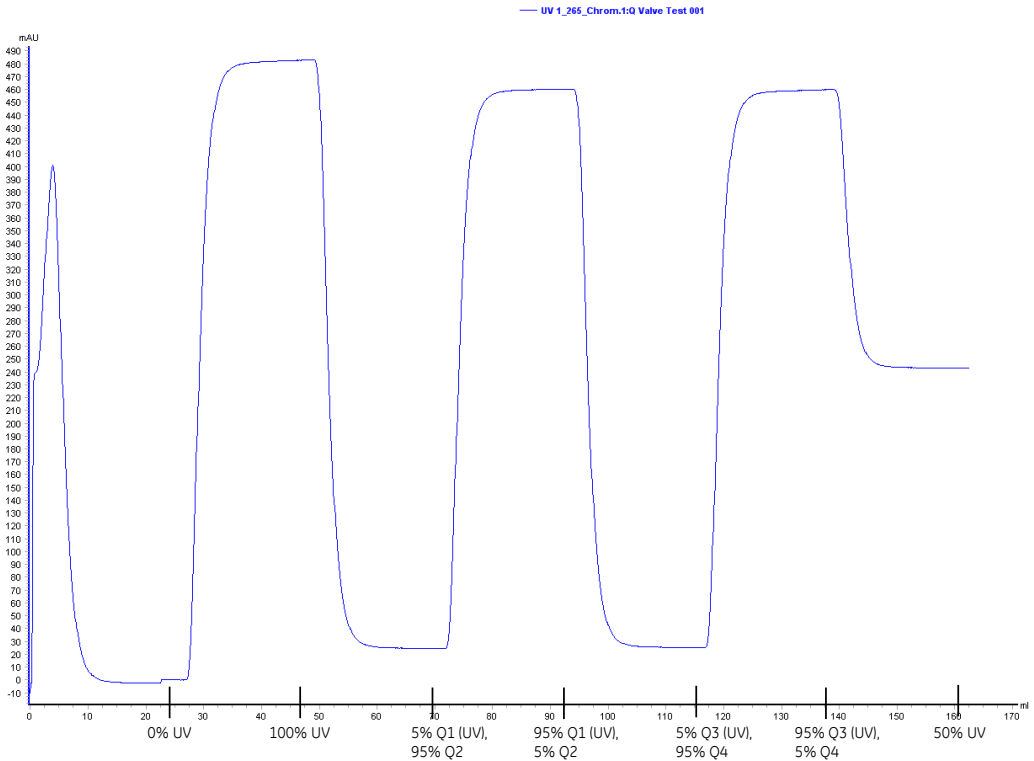

Note: The illustration shows a chromatogram from a Q valve test on ÄKTA avant 25. A chromatogram from a Q valve test on ÄKTA avant 150 has a similar appearance, but a different scale on the y-axis.

Manual evaluation - Evaluate the step response

Follow the instruction below to calculate the relative absorption plateau heights for curve UV1_265. Perform the calculations for each of the three tables in *Section 9.5 Test protocol for Q valve test, on page 107*.

Step	Action
------	--------

- | | |
|---|---|
| 1 | Subtract the baseline value (0% UV) from each of the values in column 2. Enter the results in column 3. |
|---|---|

Step	Action
2	Divide each value in column 3 by the baseline corrected value corresponding to 100% UV, multiply by 100, and enter the results in column 4.
3	Check that all values in column 4 fall within the intervals given in column 5.

Possible causes of a failed test

The table below describes possible causes of a failed test. When possible sources of error have been checked and taken care of, run the test once again.

Cause	Action
Disturbances caused by air trapped in any of the System pumps	Make sure to prime all Q inlets and purge the System pumps before method start, see <i>Section 8.3 Prime Q inlets and purge System pump A and System pump B, on page 76.</i>
Disturbances caused by damaged pump piston seals.	Replace piston seals. See <i>User Manual, chapter Maintenance.</i>
Unstable or incorrect UV signal, or drifting base line - Faulty UV monitor	See <i>User Manual, chapter Troubleshooting.</i>
Wrong mixer chamber size or faulty mixer	Replace the mixer. See <i>User Manual, chapter Maintenance.</i>
Faulty Quaternary valve	See <i>User Manual, chapter Troubleshooting.</i>
Faulty System pumps	See <i>User Manual, chapter Troubleshooting.</i>
Incorrect preparation of buffer and tubing	Make sure that the system was correctly prepared, see <i>Prepare the test, on page 90.</i>

9.4 System test

About this section

The System test checks the solvent delivery, the functionality of the UV and conductivity monitoring systems, and the valve functionality. The method run takes approximately 30 minutes. This section describes how to prepare, run, and evaluate a System test.

In this section

This section contains the following subsections:

Section	See page
9.4.1 Prepare the test	97
9.4.2 Run the test	98
9.4.3 Evaluate the test	100
9.4.4 Possible causes of a failed test	105

9.4.1 Prepare the test

Required material

The following material is required:

- Distilled water
- 1% acetone and 1 M NaCl in distilled water
- Capillary loop **Ref 1**
- Syringe, 25-30 ml
- For ÄKTA avant 25: Mixer, 1.4 ml
- For ÄKTA avant 150: Mixer, 5 ml

Prepare the test

Follow the instructions below to prepare the system before method start.

Step	Action
1	Immerse the piece of inlet tubing marked A1 in distilled water.
2	Immerse the piece of inlet tubing marked B1 in 1% acetone and 1 M NaCl in distilled water.
3	Immerse the piece of inlet tubing marked Buff in 1% acetone and 1 M NaCl in distilled water.
4	Connect the capillary loop marked Ref 1 between Column valve ports 1A and 1B to generate a backpressure.
5	For ÄKTA avant 25: Make sure that the Mixer with a chamber volume of 1.4 ml is installed. For ÄKTA avant 150: Make sure that the Mixer with a chamber volume of 5 ml is installed. For further information, see <i>User Manual chapter Maintenance</i> .
6	Prime the buffer inlets and purge System pump A and System pump B. See <i>Section 8.1 Prime buffer inlets and purge System pump A and System pump B, on page 63</i> and <i>Section 8.2 Prime sample inlets and purge the Sample pump, on page 71</i> .

9.4.2 Run the test

Start the test

Follow the instruction below to start the System test.

Step	Action
1	In the System Control module, select System:Performance Test and Report . <i>Result:</i> The System Performance Test and Report dialog opens.
2	In the System Performance Test and Report dialog. <ul style="list-style-type: none"> Select System Test Click the Run Performance Method button. <i>Result:</i> The Start Protocol of the System test opens.
3	Click Next in the dialogs of the Start Protocol to open the next dialog. The dialogs are described in table <i>Overview of the Start Protocol</i> , on page 98.
4	In the last dialog of the Start Protocol , Result name and location , click Start . <i>Result:</i> The System test starts.

Overview of the Start Protocol

The table below describes the dialogs of the **Start Protocol**.


Dialog	Description
Notes	Displays the Method Notes of the method. The Method Notes contains a method description and instructions on how to run the method. This dialog also allows the user to enter Start Notes .
Evaluation Procedures	Allows the user to select to save the report to file (recommended) and/or to print the report.
Result Name and Location	Allows the user to change result name and result location.

During the run

A **Message** dialog opens during the run. Read the messages in the dialog, and make sure that necessary preparations have been performed. Click the **Confirm and Continue** button in the **Message** dialog to change system state from **System Pause** to **Run** and proceed with the test. Alternatively, click the **Confirm** button in the **Message** dialog and click the **Continue** button on the Instrument display.

Illustration of chromatogram

The illustration below shows a chromatogram from a System test.

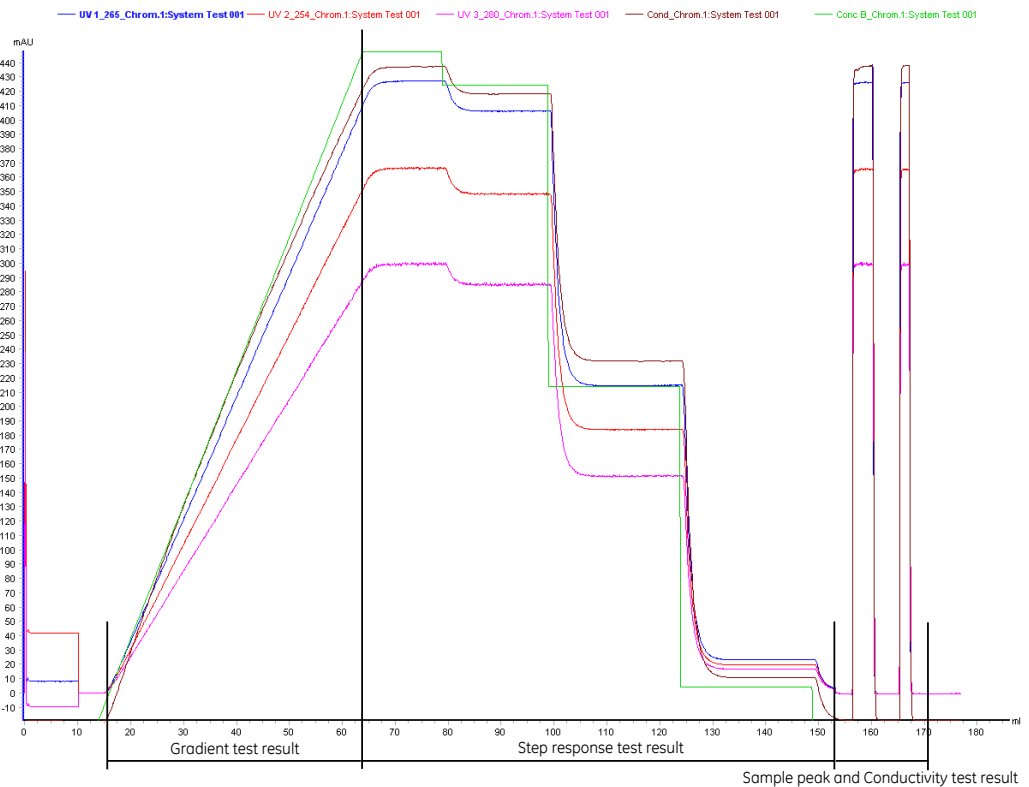

Note: The illustration shows a chromatogram from a System test on ÄKTA avant 25. A chromatogram from a System test on ÄKTA avant 150 has a similar appearance, but a different scale on the y-axis.

9.4.3 Evaluate the test

Automatic evaluation

The system automatically generates a report when the test is finished. The report can be printed in two ways:

- If the option **Save the report to file** was selected in the **Evaluation Procedures** dialog in the **Start Protocol** of the test, the report is saved in location `C:\Program Files\GE Healthcare\UNICORN\Temp`. Open the report and print it using your default Windows printer. (Recommended)

Note: *The exact search path depends on the location of the UNICORN installation folder.*

- If the option **Print report** was selected in the **Evaluation Procedures** dialog in the **Start Protocol** of the test, the report is automatically printed on the system printer. Note that a system printer must have been set up, see *UNICORN Administration and Technical Manual*.

Print the report and check the status of the tests. For each of the tests the report states "The test passed" or "The test failed".

Manual evaluation - Preparation

If you prefer to evaluate your results manually, follow the instructions below.

Step	Action
------	--------

- | | |
|---|--|
| 1 | In the Evaluation module in UNICORN, click the Open Result Navigator icon. |
|---|--|



Result: The **Result Navigator** opens.

- | | |
|---|---|
| 2 | In the Result Navigator , double-click the test result name to open the test result. |
|---|---|

- | | |
|---|----------------------------------|
| 3 | Click the Customize icon. |
|---|----------------------------------|



Result: The **Customize** dialog opens.

- | | |
|---|---|
| 4 | In the Customize dialog, select the Curves tab. |
|---|---|

Step	Action
5	<p>In the Curves tab,</p> <ul style="list-style-type: none">• Select the following curves to be displayed:<ul style="list-style-type: none">- UV1_265- UV2_254- UV3_280- Cond- Conc B• Click OK.
6	<p>In the Chromatogram pane of the Evaluation module,</p> <ul style="list-style-type: none">• Select the curve UV1_265.• Right-click and select Vertical marker.
7	<p>Click the vertical marker and drag it to the constant section of each plateau of the curve UV1_265. Read the absorbance values of each plateau in the upper left corner of the Chromatogram pane, and enter the values in the column 2 in table <i>Step response test result, on page 109</i>.</p>
8	<p>Read the absorbance values for the plateaus corresponding to 0% B and 100% B for the curves listed below. Click the curve name in Chromatogram pane to change the curve reading.</p> <ul style="list-style-type: none">• UV1_265• UV2_254• UV3_280 <p>Enter the absorbance values in column 2 in table <i>UV response test result, on page 109</i>.</p>

9 Performance tests

9.4 System test


9.4.3 Evaluate the test

Step	Action
------	--------

9	Click the Print icon.
---	------------------------------



Result: The **Print Chromatograms** dialog opens.

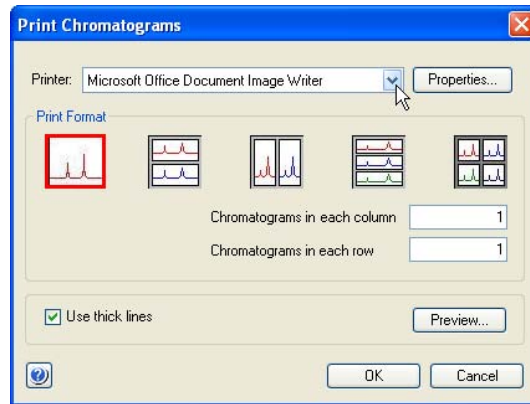

10	In the Print Chromatograms dialog:
----	---

- Select **Printer**.
- Click **OK**.

Result: Your chromatogram is printed.

Manual evaluation - Evaluate the gradient

Follow the instruction below to evaluate the gradient for curve UV1_265.

Step	Action
------	--------

- | | |
|---|--|
| 1 | Place a ruler along the gradient part of curve UV1_265 in the printed report. |
| 2 | Check that the curve is linear between 5% B and 90% B and void of discontinuities. |
| 3 | Enter the start and stop values of the linear section of the curve in <i>Gradient test result</i> , on page 109. |

Manual evaluation - Evaluate the step response

Follow the instruction below to calculate the relative absorption plateau heights for curve UV1_265.

Step	Action
1	Subtract the baseline value (0% B) from each of the values in column 2 in table <i>Step response test result</i> , on page 109. Enter the results in column 3.
2	Divide each value in column 3 by the baseline corrected value corresponding to 100% B, multiply by 100, and enter the results in column 4.
3	Check that all values in column 4 fall within the intervals given in column 5.

Manual evaluation - Evaluate the UV response

Follow the instruction below to calculate the UV response.

Step	Action
1	Subtract the baseline values (0% B) corresponding to each UV curve from the values corresponding to 100% B in column 2 in table <i>UV response test result</i> , on page 109. Enter the results in column 3.
2	Calculate the absorbance ratios 265/254 nm and 265/280 nm using the values in column 3. Enter the values in column 4.
3	Check that the obtained ratios fall within the intervals given in column 5.

Manual evaluation - Evaluate the Sample Peak and Conductivity Test Result

Follow the instruction below to inspect the two sample peaks of the Cond curve in the chromatogram.

Step	Action
1	Check that the Cond curve has a shape similar to the shape of the UV curves.

9 Performance tests

9.4 System test

9.4.3 Evaluate the test

Step	Action
2	<p>Check that the baseline of the Cond curve looks correct and is not drifting.</p> <ul style="list-style-type: none">• If the baseline is not drifting: Integrate the curve using the function Zero baseline.• If the baseline is drifting: Correct the baseline before integrating the curve. <p>For further information on peak integration, see <i>Evaluation Manual</i>.</p>
3	<p>Enter the area of each peak in column 2 in table <i>Sample peak and conductivity test result</i>, on page 110. Add the areas of peak 1 and peak 2, and divide the area of each peak with the total peak area. Enter the results in column 3.</p>
4	<p>Check that the obtained percentages fall within the intervals given in column 4.</p>

9.4.4 Possible causes of a failed test

Introduction

The tables in this section describe possible causes of a failed test. When possible sources of error have been checked and taken care of, run the test once again.

Faulty Gradient Test Result

Cause	Action
Disturbances caused by air trapped in any of the pumps	Make sure to prime the buffer inlets and to purge the System pumps and the Sample pump before method start. See <i>Section 8.1 Prime buffer inlets and purge System pump A and System pump B, on page 63</i> and <i>Section 8.2 Prime sample inlets and purge the Sample pump, on page 71</i> .
Disturbances caused by damaged pump piston seals.	Replace piston seals. See <i>User Manual chapter Maintenance</i> .
Unstable or incorrect UV signal, or drifting base line - Faulty UV monitor	See <i>User Manual chapter Troubleshooting</i>
Wrong Mixer chamber size or faulty Mixer	Replace the Mixer. See <i>User Manual chapter Maintenance</i> .

Faulty Step Response Result

Cause	Action
If all values are faulty - air in the pump or a faulty pump	Air in pumps: Make sure to prime the buffer inlets and to purge the System pumps and the Sample pump before method start. See <i>Section 8.1 Prime buffer inlets and purge System pump A and System pump B, on page 63</i> and <i>Section 8.2 Prime sample inlets and purge the Sample pump, on page 71</i> . Faulty pump: See <i>User Manual chapter Troubleshooting</i> .
Faulty values at 5% - damaged pump piston seal in System pump B	Replace piston seals. See <i>User Manual chapter Maintenance</i> .

9 Performance tests

9.4 System test

9.4.4 Possible causes of a failed test

Cause	Action
Faulty values at 95% - damaged pump piston seal in System pump A	Replace piston seals. See <i>User Manual chapter Maintenance</i> .

Faulty UV Response Test Result

Cause	Action
Faulty UV monitor	Restart the instrument to calibrate the UV monitor. See <i>User Manual chapter Troubleshooting</i> .

Faulty Runlog Event Test Result

Cause	Action
Any of the valves do not switch	See <i>User Manual chapter Troubleshooting</i> .

Faulty Sample Peak and Conductivity Test Result

Cause	Action
If the shape of the Cond curve differs from the shape of the UV curves - faulty Conductivity monitor	See <i>User Manual, chapter Troubleshooting</i> .
Unproportional peaks - faulty Sample pump	See <i>User Manual, chapter Troubleshooting</i> .
Faulty baseline	Perform a manual peak integration of the Cond curve. See <i>Evaluation Manual</i> .

9.5 Test protocol for Q valve test

Step response test result - Q1, Q2

1 Programmed Conc.	2 Value read (mAU)	3 Baseline corrected value	4 Normalized value	5 Allowed interval
100% UV (100% Q1)				
95% UV (95% Q1, 5% Q2)				93.6 - 96.4
5% UV (5% Q1, 95% Q2)				3.6 - 6.4
0% UV (100% Q2)				

Step response test result - Q3, Q4

1 Programmed Conc.	2 Value read (mAU)	3 Baseline corrected value	4 Normalized value	5 Allowed interval
100% UV (100% Q1)				
95% UV (95% Q3, 5% Q4)				93.6 - 96.4
5% UV (5% Q3, 95% Q4)				3.6 - 6.4
0% UV (100% Q2)				

Step response test result

- Q1, Q2, Q3, Q4

1 Programmed Conc.	2 Value read (mAU)	3 Baseline corrected value	4 Normalized value	5 Allowed interval
100% UV (100% Q1)				
50% UV (25% Q1, 25% Q2, 25% Q3, 25% Q4)				48.6 - 51.4
0% UV (100% Q2)				

9.6 Test protocol for System test

Gradient test result

Gradient linear from % B to % B. (The gradient must be linear from 5% to 90% with an allowed deviation of 1.00%.)

Step response test result

1 Programmed Conc. %B	2 Value read (mAU)	3 Baseline corrected value	4 Normalized value	5 Allowed interval
100			-	-
95				94.00 - 96.00
50				49.00 - 51.00
5				4.00 - 6.00
0			-	-

UV response test result

1 Wavelength (nm)	2 Value read (mAU)		3 Baseline corrected value	4 Normalized value	5 Allowed interval
	100% B	0% B			
254				-	-
265/254	-	-	-		1.10 - 1.23
265				-	-
265/280	-	-	-		1.30 - 1.58
280				-	-

Sample peak and conductivity test result

1 Peak no.	2 Area	3 Area/ Peak area (%)	4 Allowed interval (%)
1			65.00-67.00
2			33.00-35.00

10 Reference information

About this chapter

This chapter lists the technical specifications of the ÄKTA avant. The chapter also includes a list of wetted materials and a chemical resistance guide.

In this chapter

This chapter contains the following sections:

Section	See page
10.1 System specifications	112
10.2 Component specifications	114
10.3 Wetted materials	121
10.4 Chemical resistance guide	124

10.1 System specifications

Introduction

This section lists the system specification data of ÄKTA avant. For components data see *Section 10.2 Component specifications, on page 114.*

System specifications

Parameter	Data
System configuration	Benchtop system, external computer
Control system	For ÄKTA avant 25: UNICORN 6.0 or higher version For ÄKTA avant 150: UNICORN 6.1 or higher version
Connection between PC and instrument	Ethernet
Dimensions (W x D x H)	860 x 710 x 660 mm
Weight (excluding computer)	116 kg
Power supply	100-240 V ~, 50-60 Hz
Power consumption	800 VA
Enclosure protective class	IP 21, wet side IP 22

Parameter	Data
Tubing and connectors	<p>ÄKTA avant 25:</p> <ul style="list-style-type: none"> • Inlet: FEP tubing, i.d. 1.6 mm, Tubing connector 5/16" + Ferrule (yellow), 1/8" • Pump to Injection valve: PEEK tubing, i.d. 0.75 mm, Fingertight connector, 1/16" • After Injection valve: PEEK tubing, i.d. 0.50 mm, Fingertight connector, 1/16" <p>ÄKTA avant 150:</p> <ul style="list-style-type: none"> • Inlet: FEP tubing, i.d. 2.9 mm, Tubing connector 5/16" + Ferrule (blue), 3/16" • After pumps: PEEK tubing, 1.0 mm i.d., Fingertight connector, 1/16"

Environmental ranges

Parameter	Data
Storage and transport temperature range	-25°C to +60°C
Chemical environment	See <i>Section 10.4 Chemical resistance guide</i> , on page 124.

Operating range

Parameter	Data
Operating temperature range	4°C to 35°C
Relative humidity	20% to 95%, non-condensing

10.2 Component specifications

Introduction

This section specifies the operating data of the components in ÄKTA avant. For general data for the system see *Section 10.1 System specifications, on page 112*.

System pumps

Parameter	Data
Pump type	Piston pump, metering type
Flow rate range	ÄKTA avant 25: 0.001 to 25 ml/min ÄKTA avant 150: 0.01 to 150 ml/min
Pressure range	ÄKTA avant 25: 0 to 20 MPa (2900 psi) ÄKTA avant 150: 0 to 5 MPa (725 psi)
Viscosity range	ÄKTA avant 25: 0.35 to 10 cP ÄKTA avant 150: 0.35 to 5 cP
Flow rate specifications	ÄKTA avant 25: <ul style="list-style-type: none">• Accuracy: $\pm 1.2\%$• Precision: RSD < 0.5% (Conditions: 0.25 to 25 ml/min, < 3 MPa, 0.8 to 2 cP) ÄKTA avant 150: <ul style="list-style-type: none">• Conditions: 1.0 to 150 ml/min, < 3 MPa, 0.8 to 2 cP• Accuracy: $\pm 1.5\%$• RSD < 0.5%

Sample pump

Parameter	Data
Pump type	Piston pump, metering type
Flow rate range	ÄKTA avant 25: 0.01 to 25 ml/min ÄKTA avant 150: 0.01 to 150 ml/min
Pressure range	ÄKTA avant 25: 0 to 10 MPa (1450 psi) ÄKTA avant 150: 0 to 5 MPa (725 psi)
Viscosity range	0.7 to 10 cP
Flow rate specifications	<p>ÄKTA avant 25:</p> <ul style="list-style-type: none"> • Conditions: 0.25 to 25 ml/min, < 3 MPa, 0.8 to 3 cP • Accuracy: $\pm 2\%$ • RSD < 0.5% <p>ÄKTA avant 150:</p> <ul style="list-style-type: none"> • Conditions: 1.0 to 150 ml/min, < 3 MPa, 0.8 to 2 cP • Accuracy: $\pm 2\%$ • RSD < 0.5%

Valves

Parameter	Data
Type	Rotary valves
Number of valves	Up to 12
Functions	Inlet A, Inlet B, Sample inlet, Injection, Column, pH, Outlet

Quaternary valve

Parameter	Data
Type	4-port solenoid actuated membrane valve
Functions	Quaternary gradients or BufferPro

Number of inlets

Parameter	Data
Inlet A	7, expandable to 14
Inlet B	7, expandable to 14
Sample inlet	7, expandable to 14
Quaternary inlet	4, expandable to 18

Pressure sensors

Parameter	Data
Placement of sensors	System pump, Sample pump, Pre-Column, Post-Column
Range	ÄKTA avant 25: 0 to 20 MPa (2900 psi) ÄKTA avant 150: 0 to 5 MPa (725 psi)
Accuracy	ÄKTA avant 25: ± 0.02 MPa or $\pm 2\%$ whichever is greater ÄKTA avant 150: ± 0.015 MPa or $\pm 1.5\%$ whichever is greater

Air sensors

Parameter	Data
Placement of sensors	Inlet A, Inlet B, Sample inlet

Parameter	Data
Optional placement	ÄKTA avant 25: Before Sample inlet valve or after Injection valve ÄKTA avant 150: After Injection valve
Sensing principle	Ultrasonic

Outlet valve fractionation

Parameter	Data
Number of outlets	10, expandable to 32
Fraction volumes	ÄKTA avant 25: 0.01 to 100 000 ml ÄKTA avant 150: 1 to 100 000 ml
Delay volume (UV – Outlet valve)	ÄKTA avant 25: 142 µl ÄKTA avant 150: 535 µl

Mixer

Parameter	Data
Mixing principle	Chamber with magnetic stirrer
Mixer volume	ÄKTA avant 25: 0.6, 1.4 or 5 ml ÄKTA avant 150: 1.4, 5 or 15 ml

Gradient formation

Parameter	Data
Gradient flow rate range	ÄKTA avant 25: <ul style="list-style-type: none"> • Binary: 0.25 to 25 ml/min • Quaternary: 0.5 to 25 ml/min ÄKTA avant 150: <ul style="list-style-type: none"> • Binary: 1 to 150 ml/min • Quaternary: 2 to 40 ml/min
Gradient composition accuracy	Binary: $\pm 0.6\%$ Quaternary: $\pm 1\%$

Pressure monitors

Parameter	Data
Placement of sensors	System pump
Range	0 to 20 MPa (2900 psi)
Accuracy	± 0.02 MPa or $\pm 2\%$ whichever is greater

UV monitor

Parameter	Data
Wavelength range	190 to 700 nm in steps of 1 nm, up to 3 wavelengths
Absorbance range	-6 to +6 AU
Linearity	within $\pm 2\%$ at 0 to 2 AU
Operating pressure	0 to 2 MPa (290 psi)
Flow cells	<ul style="list-style-type: none"> • 0.5 mm optical path length, 1 μl cell volume • 2 mm optical path length, 2 μl cell volume • 10 mm optical path length, 8 μl cell volume

Conductivity monitor

Parameter	Data
Conductivity reading range	0.01 mS/cm to 999.99 mS/cm
Accuracy	± 0.01 mS/cm or ± 2%, whichever is greater, (within 0.3 to 300 mS/cm)
Operating pressure	0 to 5 MPa (725 psi)
Flow cell volume	22 µl

Temperature monitor

Parameter	Data
Reading range	0°C to 99°C
Accuracy	± 1.5°C within 4°C to 45°C

pH monitor

Parameter	Data
pH reading range	0 to 14
Accuracy	± 0.1 pH unit (within pH 2 to 12, temp. within 3°C from calibration temp.)
Operating pressure	0 to 0.5 MPa (72 psi)
Flow cell volume	76 µl

Fraction collector

Parameter	Data
Number of fractions	up to 576

10 Reference information

10.2 Component specifications

Parameter	Data
Vessel types	3, 8, 15 or 50 ml tubes 250 ml bottles Deep well plates, 96 / 48 / 24
Vessel type selection	Automatic recognition
Fraction volumes	0.1 to 50 ml
Spillage-free modes	ÄKTA avant 25: Drop sync or Accumulator ÄKTA avant 150: Accumulator
Protection of fractions	Covered vessels and temperature control
Organic solvents	Ei
Delay volume (UV – Dispenser head)	ÄKTA avant 25: 518 µl ÄKTA avant 150: 1831 µl
Temperature control	Target temperature: 6°C to 20°C

10.3 Wetted materials

Introduction

This section specifies the wetted materials of ÄKTA avant.

Wetted materials

The following wetted materials are used in ÄKTA avant:

Material	Component
Aluminium oxide	Sample pump
Borosilicate	Fraction collector
ECTFE, EthyleneChloroTriFluoroEthylene	System pumps Sample pump
ETFE, EthyleneTetraFluoroEthylene	Tubing
EPDM, EthylenePropyleneDiene M-class rubber	Quaternary valve Fraction collector
Elgiloy™/UHMWPE, UltraHighMolecularWeightPolyEthylene	System pumps Sample pump
FEP, FluorinatedEthylenePropylene	Tubing
FFPM/FFKM, FullyFluorinatedPropyleneMonomer	pH valve Flow restrictor Pressure monitors
FPM/FKM, FluorinatedPropyleneMonomer	Mixer chamber Sample pump
Hastelloy™ C-276	Pressure monitors
PCTFE, PolyChloroTriFluoroEthylene	Conductivity monitor

10 Reference information

10.3 Wetted materials

Material	Component
PEEK, PolyEtherEtherKetone	Quaternary valve Inlet valves System pumps Pressure monitors Mixer chamber Injection valve Column valve UV flow cell Conductivity monitor pH valve Flow restrictor Sample pump Outlet valve Fraction collector Tubing
PolyEthylene	Fraction collector
PolyImide	UV flow cell
PP, PolyPropylen	Mixer chamber Fraction collector
PPS, PolyPhenyleneSulfide	Fraction collector
PTFE, PolyTetraFluorEthylene	Mixer chamber Injection valve Column valve UV flow cell pH valve Outlet valve
PVDF, PolyVinylideneDiFluoride	System pumps Sample pump
Ruby	System pumps Sample pump

Material	Component
Sapphire	System pumps Sample pump
Silica	UV flow cell
Silicone	Tubing
Titanium	System pumps Pressure monitors Conductivity monitor Sample pump

10.4 Chemical resistance guide

Introduction

This section specifies the chemical resistance of ÄKTA avant to some of the most commonly used chemicals in liquid chromatography.

Biocompatibility

ÄKTA avant is designed for maximum biocompatibility, with biochemically inert flow paths constructed mainly from titanium, PEEK and highly resistant fluoropolymers and fluoroelastomers. Titanium is used as far as possible to minimize contribution of potentially deactivating metal ions such as iron, nickel and chromium. There is no standard stainless steel in the flow path. Plastics and rubber materials are selected to avoid leakage of monomers, plasticizers or other additives.

Cleaning chemicals overview

Strong cleaning works well with 2 M sodium hydroxide, 70% acetic acid or the alcohols methanol, ethanol and isopropyl alcohol. Complete system cleaning with 1 M hydrochloric acid should be avoided due to sensitivity in the pressure sensors. For cleaning separation media with 1 M hydrochloric acid, use loop injections of the acid. Make sure that the column is not mounted on a column valve **V9-C** (which contains a pressure sensor).

If sodium hypochlorite is used as sanitizing agent instead of 2 M sodium hydroxide, use a concentration up to 10%.

Chromatography chemicals overview

Reversed phase chromatography of proteins works well with 100% acetonitrile and additives trifluoroacetic acid (TFA) up to 0.2% or formic acid up to 5%.

Strong organic solvents like ethyl acetate, 100% acetone or chlorinated organic solvents should be avoided. These might cause swelling of plastic material and reduce the pressure tolerance of PEEK tubing. For this reason, flash chromatography and straight ("normal") phase is generally not recommended on the system

Assumptions made

The ratings are based on the following assumptions:

- Synergy effects of chemical mixtures have not been taken into account.
- Room temperature and limited overpressure is assumed.

Note: *Chemical influences are time and pressure dependent. Unless otherwise stated, all concentrations are 100%.*

List of chemicals

Proposed chemical compatibility for ÄKTA avant. All chemicals used for CIP and cleaning are for short term use only, ambient temperature < 25°C, if not other stated.

Note: *A user can be exposed to large volumes of chemical substances over a long time period. A Material Safety Data Sheet (MSDS) provides the user with information regarding characteristics, human and environmental risks and preventive measures. Make sure that you have the MSDS available from your chemical distributor and/or databases on the internet.*

Chemical	Concentration	CAS no/ EC no	Usage
Aqueous buffers, pH 2-12	-	N/A	Separation
Acetic acid	70%	64-19-7/ 200-580-7	Cleaning In Place (CIP)
Acetonitrile ¹ Note: <i>Quaternary valve is not resistant. Depending on pressure, tubing between pump head and pressure monitor needs to be changed. See User Manual, chapter Prepare the system for a run.</i>	100%	75-05-8/ 200-835-2	Reversed Phase Chromatography (RPC)
Acetonitrile/THF Note: <i>Quaternary valve is not resistant.</i>	85/15	109-99-9/ 203-726-8 (Tetrahydrofuran)	RPC

10 Reference information
10.4 Chemical resistance guide

Chemical	Concentration	CAS no/ EC no	Usage
Acetone	10%	67-64-1/ 200-662-2	Rare: RPC
Ammonia	30%	7664-41-7/ 231-635-3	Oligonucleotide synthesis
Ammonium chloride	2 M	12125-02-9/ 235-186-4	Hydrophobic Interaction Chromatography (HIC)
Ammonium sulphate	3 M	7783-20-2/ 231-984-1	Cleaning of plasmids, HIC
Arginine	2 M	74-79-3/ 200-811-1	Wash and elution (pH 4), using protein A media, refolding
Benzyl alcohol	2%	100-51-6/ 202-859-9	Cleaning and storage of columns
Decon™ 90	10%	N/A	Cleaning
Dimethyl sulphoxide (DMSO)	5%	67-68-5/ 200-664-3	CIP, RPC, cell separation
Dithiothreitol (DTT)	100 mM	3483-12-3 / 222-468-7	Reducing agent
Dithioerythritol (DTE)	100 mM	6892-68-8/ 229-998-8	Reducing agent
Ethylenediaminetetraacetic acid (EDTA)	100 mM	60-00-4/ 200-449-4	Buffer additive
Ethanol	20%	75-08-1/ 200-837-3	Storage Note: Long term use (1 month)
Ethanol	96%	75-08-1/ 200-837-3	CIP
Ethanol + NaOH	40% + 1 M	N/A	CIP
Ethylene glycol	50%	107-21-1/ 203-473-3	HIC elution
Formic acid	1%	64-18-6/ 200-579-1	RPC, peptide separation
Glycerol	50%	56-81-5/ 200-289-5	HIC

Chemical	Concentration	CAS no/ EC no	Usage
Glycine	2 M	56-40-6/ 200-272-2	Cleaning of MAb binding media
Guanidinium hydrochloride	6 M	50-01-1/ 200-002-3	Denaturizing agent
Hydrochloric acid ²	max. 0.1 M	7647-01-0/ 231-595-7	CIP
Imidazole	2 M	288-32-4/ 206-019-2	2 M needed for buffer preparation
Isopropanol	100%	67-63-0/ 200-661-7	CIP fatty gels
Methanol	100%	74-93-1/ 200-659-6	RPC, CIP
Mercaptoethanol	20 mM	37482-11-4/ 253-523-3	Reducing agent
n-Propanol	50%	67-63-0/ 200-661-7	RPC
Phosphoric acid	0.1 M	7664-38-2/ 231-633-2	CIP
Potassium phosphate	1 M	7778-77-0/ 231-913-4	HIC
Potassium chloride	4 M	7447-40-7/ 231-211-8	4M for pH cell storage Note: Long term use (1 month)
Sodium dodecyl sulfate (SDS)	1%	151-21-3/ 205-788-1	Detergent
Sodium chloride	4 M	7647-14-5/ 231-598-3	CIP
Sodium hydroxide	2 M	1310-73-2/ 215-185-5	CIP
Sodium sulphate	1 M	7757-82-6/ 231-820-9	HIC
Trichloroacetic acid	1%	76-03-9/ 200-927-2	RPC additive
Trifluoroacetic acid	1%	176-05-1/ 200-929-3	RPC additive, peptide separation
Triton™-X	1%	9002-93-1 (Triton X 100)	Detergent
Tween™	1%	9005-64-5/ 500-018-3 (Tween 20)	Detergent

Chemical	Concentration	CAS no/ EC no	Usage
Urea	8 M	57-13-6/ 200-315-5	Buffer additive, denaturing agent
Water	100%	7732-18-5/ 231-791-2	Note: Long term use (1 month)

- 1 PEEK tubing is biocompatible and chemically inert to most solvents used in purification of proteins. It has in general very good pressure limits, especially for water based buffers. However, organic solvents penetrates weaknesses in the tubing walls easier than water based buffers. Special care should therefore be taken with prolonged use of organic solvents close to pressure limits.
- 2 If hydrochloric acid, HCl, is used as a cleaning agent when columns are connected to the system, the HCl concentration should not exceed 0.1 M in the pressure sensors. Remember that the ÄKTA avant system has pressure sensors in the Column valve V9-C.
For other parts of the system up to 1 M HCl is acceptable for short time use.
Long time use of 0.2 M HCl connected to the Quaternary valve as part of a **BufferPro** recipe is acceptable. The solution becomes diluted further down in the system.

For local office contact information, visit
www.gelifesciences.com/contact

GE Healthcare Bio-Sciences AB
Björkgatan 30
751 84 Uppsala
Sweden

www.gelifesciences.com/AKTA

GE, imagination at work and GE monogram are trademarks of General Electric Company.

ÄKTA, UNICORN, and Whatman are trademarks of GE Healthcare companies.

Eppendorf is a trademark of Eppendorf AG.

Hastelloy is a trademark of Hayne's Intl, Inc.

Seahorse Bioscience is a trademark of Seahorse Bioscience Inc.

Microsoft and Windows are trademarks of Microsoft Corporation.

© 2009-2012 General Electric Company – All rights reserved.

First published Dec. 2012

All goods and services are sold subject to the terms and conditions of sale of the company within GE Healthcare which supplies them. A copy of these terms and conditions is available on request. Contact your local GE Healthcare representative for the most current information.

UNICORN: Any use of this software is subject to GE Healthcare Standard Software End-User License Agreement for Life Sciences Software Products. A copy of this Standard Software End-User License Agreement is available on request.

GE Healthcare Europe GmbH
Munzinger Strasse 5, D-79111 Freiburg, Germany

GE Healthcare UK Limited
Amersham Place, Little Chalfont, Buckinghamshire, HP7 9NA, UK

GE Healthcare Bio-Sciences Corp.
800 Centennial Avenue, P.O. Box 1327, Piscataway, NJ 08855-1327, USA

GE Healthcare Japan Corporation
Sanken Bldg. 3-25-1, Hyakunincho Shinjuku-ku, Tokyo 169-0073, Japan



imagination at work