

User Instruction

calaject™

Computer Assisted Local Analgesia



RØNVIG Dental Mfg. A/S

USER INSTRUCTION

CONGRATULATIONS ON YOUR NEW CALAJECT™!

Please read these instructions thoroughly before you start using your CALAJECT™



CALAJECT™ MAY ONLY BE USED BY TRAINED PERSONNEL authorized to perform dental injections. For this reason, this manual does not include specific instructions in injection techniques. The manufacturer cannot be held liable for injuries in patients due to unauthorized or incorrect use.


CONTENTS

- | | |
|-----------------------|-----------------------|
| 1 Control unit | 1 Stand for handpiece |
| 1 Handpiece with cord | 1 Charger |
| 1 Footswitch | 1 User instructions |
| 3 Cartridge barrels | |

DESCRIPTION OF CONTROL UNIT


FRONT PANEL WITH FINGER TOUCH DISPLAY

1.  On/Off switch
2. Bar-scale for display of current injection pressure/resistance
3. Program selection 1, 2, 3
4.  Piston rod retraction. Returns the piston rod to start position
5. Charging and battery level indicator

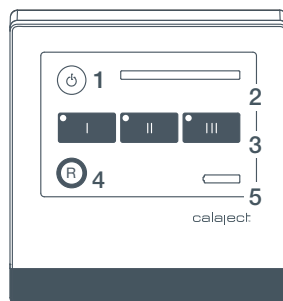

REAR PANEL

6. Charger socket
7. Socket for handpiece
8. Connection to foot control or footswitch of the dental unit. Some multi-function footswitches have a spare function (e.g. a call function) that can be used for operating CALAJECT™
9. Volume control for sound signal
10. Sound aperture

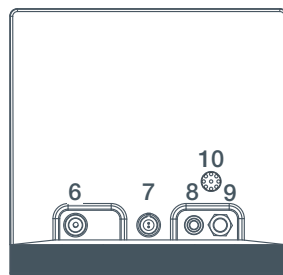

RECOMMENDATION

It is recommended that new users carry out test injections in the air in order to familiarise themselves with the three different programs.

CLEANING

- The CALAJECT™ unit, the handpiece and the handpiece stand can be cleaned with a pad moistened with any surface disinfectant used in dental practice.



Do not immerse into liquid. Do not autoclave. The CALAJECT™ unit and the handpiece contain sensitive electronic components that do not withstand sterilization or immersion in liquid.

- The cartridge barrel can be autoclaved at max. 135°C.

The cartridge barrel may get a frosted appearance after a number of sterilizations. That will not affect the strength of the cartridge barrel, but it is recommended to replace the barrel regularly in order to keep a clear view to the cartridge.

If the cartridge barrel is damaged it should be replaced. Non-original barrels may not be used.

Replacement barrels (3 pcs. per pack) can be ordered over your CALAJECT™ dealer.

SERVICE WARRANTY & REPAIR

CALAJECT™ is covered by a 2-year guarantee on materials and construction.

Normal wear and tear and damages due to inadequate use or maintenance are not covered by the warranty.

In the event of malfunction, please return the device to your CALAJECT™ dealer for repair.

OPTION

CALAJECT™ is operated by a separate foot control included in the package. But the system can also be connected to the multi-function foot switch of your dental unit by your dental service engineer. When connecting to the dental unit choose a potential free contact (relaycontact). The assistant call function can often be used for this purpose. Connect to Calaject with a screened cable less than 3m long and with a 3.5mm minijack plug.

GETTING STARTED

- **Connect the handpiece cord plug to the CALAJECT™ rear panel**
 - the red dot on CALAJECT™ and the handpiece plug must be aligned.
- Unplug by pulling the grooved sliding ring backwards (do not turn).



CALAJECT™ should not be placed close to devices that are sensitive to
- or generate - electromagnetic interference.

- **Check the battery status on the display.** Charging time approx. 3 hours.
Operating time approx. 5 hours.
- **Attach a needle on the cartridge barrel and insert an anaesthetic cartridge.**
To avoid leakage at the membrane of the cartridge it is advised to screw on the needle first and then insert the cartridge. The cartridge barrel fits standard 1,7/1,8 ml dental cartridges and standard dental needles.
- **Mount the loaded cartridge barrel onto the handpiece.** Before the barrel is screwed onto the handpiece, the piston rod should be retracted to starting position. It will automatically return to starting position, when CALAJECT™ is turned on. The piston rod can also be returned to starting position by pressing "R" on the display.
- **Activate the foot control,** until the anaesthetic solution is seen to come out of the needle.
- **Select program.** When the foot control is activated again, the chosen injection program will be active.
- CALAJECT™ will stop automatically, when the cartridge is empty.
Return the piston rod by pressing "R" at the display.


CALAJECT™ will stop automatically when it reaches the pre-programmed maximum pressure. In such case, a long sound signal will be heard and the pressure bar-scale at the display will turn off. Wait a moment or move the needle to a new position before you continue the injection.

PROGRAM 1

Recommended for intraligamental – and also palatal – analgesia

- Activate the foot control → slow injection speed, approx. 0,006 ml/sec.
Optional: Release / re-activate the foot control → the injection speed will increase to 0,009 ml/sec.

The PDL-technique requires a relatively high injection pressure initially. This is why Program 1 allows a substantially higher injection pressure / resistance than program 2 and 3 before the automatic safety stop is activated.

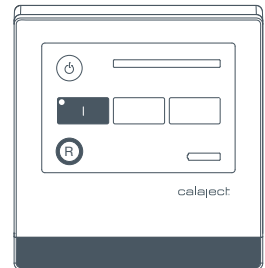

TIP For intraligamental (PDL) analgesia it is recommended to dose 0,2 – 0,9 ml per root depending on the size of the root and the expected duration of the procedure. For further guidance on the PDL technique, we refer to the published literature on the subject.


TIP If the pressure has become so high that CALAJECT™ stops, the needle opening may have been blocked and it is recommended to rotate the needle slightly in order to obtain a good flow.

TIP Auto pilot - five seconds after program start the sound signal is changing. This indicates that you can release the foot control and let the auto pilot take over. You interrupt the injection by re-activating the foot control. Note, auto pilot is only an option in program 1.

PROGRAM 2

Recommended for infiltration analgesia

- Initially 10 seconds with slow injection speed (approx. 0,006 ml/sec). During the subsequent 5 seconds it will gradually increase to medium injection speed of 0,03 ml/sec.
- Aspirates automatically whenever the foot control is released. The small back-suction will also prevent after-dripping from the needle.

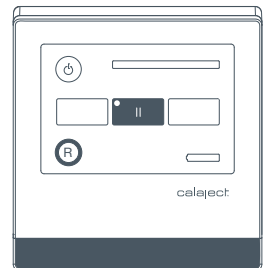

PROGRAM 3

Recommended for regional nerve block analgesia

- Initially slow injection speed (approx. 0,006 ml/sec). By releasing/re-activating the foot control in one swift movement, the injection speed will increase gradually over the next 5 seconds to high speed (approx. 0,04 ml/sec). Hereafter, high speed at every stop/start.
- Aspirates automatically whenever the foot control is released.

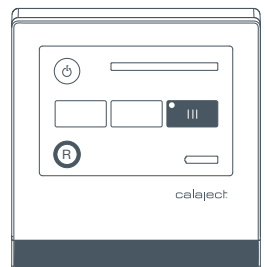


Table 1 Electromagnetic emissions		
The “CALAJECT” is intended for use in the electromagnetic environment specified below. The customer or the user of the “CALAJECT” should assure that it is used in such an environment.		
Emissions test	Compliance	Electromagnetic environment guidance
RF emissions CISPR11	Group 1	The device uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR11	Class B	The device is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
Harmonic emissions IEC 61000-3-2	Class A	
Voltage fluctuations/ flicker emissions IEC 61000-3-3	Complies	

Table 2 Electromagnetic immunity			
The “CALAJECT” is intended for use in the electromagnetic environment specified below. The customer or the user of the “CALAJECT” should assure that it is used in such an environment.			
Immunity test	IEC 60601 Test level	Compliance level	Electromagnetic environment guid- ance
Electrostatic discharge (ESD) IEC61000-4-2	±6KV contact ±8KV air	±6KV contact ±8KV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 40%.
Electrical fast transient/ burst IEC61000-4-4	±2KV for power supply lines ±1KV for input/output lines	±2KV for power supply lines ±1KV for input/output lines	Mains power supply quality should be that of typical residential area.
Surge IEC61000-4-5	±1KV differential mode ±2KV common mode	±1KV differential mode ±2KV common mode	
Voltage dips, short interruptions and voltage variations on power supply input lines. IEC61000-4-11	<5% UT for 0.5 cycle 40%UT for 5 cycles 70%UT for 25 cycles <5%UT for 5 seconds	<5% UT for 0.5 cycle 40%UT for 5 cycles 70%UT for 25 cycles <5%UT for 5 seconds	Mains power supply quality should be that of typical residential area.
Power frequency (50-60Hz) magnetic field. IEC61000-4-8	3A/m	3A/m	Power frequency magnetic fields should be that of typical residential area.

Specifications:

Charge:: 10-14v/1ADC, medical
Output:: 0-4v/0,3ADC

Temperature range:

Operating: 10 - 35°C
Storage: -20 – 60°C
Humidity: 10 – 95%

Classifikations: COUNCIL DIRECTIVE 93/42/EEC Class IIa

Standards: EN60601-1

Disposal: Separate collection for electronic equipment


Table 3 Electromagnetic immunity			
The “CALAJECT” is intended for use in the electromagnetic environment specified below. The customer or the user of the “CALAJECT” should assure that it is used in such an environment.			
Immunity test	IEC 60601 Test level	Compliance level	Electromagnetic environment guidance
Conducted RF IEC61000-4-6	3Vrms 150KHz to 80 MHz	3Vrms	Portable and mobile RF communications equipment should be used no closer to any part of the device, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter. Recommended separation distance $d=1,17\sqrt{P}$ $d=1,17\sqrt{P}$ 80 MHz to Hz to 800 MHz $d=2,23\sqrt{P}$ 800 MHz to 2,5 GHz where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in meters (m). Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey ^a , should be less than the compliance level in each frequency range. ^b Interference may occur in the vicinity of equipment marked with the following symbol: 
Radiated RF IEC61000-4-3	3V/can 80MHz to 2,5GHz	3V/m	
NOTE 1 At 80 MHz and 800 MHz, the higher frequency range applies.			
NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.			
^a Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the device is used exceeds the applicable RF compliance level above, the device should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the device.			
^b Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 3 V/m.			

Table 4 Recommended separation distances between portable and mobile RF communications equipment and the “CALAJECT”			
The device is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the device can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the device as recommended below, according to the maximum output power of the communications equipment			
Rated maximum output of transmitter (W)	Separation distance according to frequency of transmitter m		
	150 kHz to 80 MHz $d = 1,17 P$	80 MHz to 800 MHz $d = 1,17 P$	800 MHz to 2,5 GHz $d = 3,5 P$
0,01	0,12	0,12	0,35
0,1	0,37	0,37	1,11
1	1,2	1,2	3,5
10	3,7	3,7	11,1
100	12	12	35
For transmitters rated at a maximum output power not listed above the recommended separation distance d in metres (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.			
NOTE 1 At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.			
NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.			

TECHNICAL SPECIFICATIONS			
	Control unit	Handpiece	Stand
Length	95 mm	200 mm (incl. barrel)	
Width	120 mm	Ø 12 mm	Ø 60 mm
Height	115 mm		34 mm
Weight	750 g	50 g	410 g
Nominal voltage	90-240 V – 50/60Hz		
Battery (Lithium-ion)	8 hours on each charge		
Charging time	Approx. 3 hours		
Dental Cartridge		1,7/1,8 ml standard cartridges	
Dental Needles		Standard M6 and 7/32”	

Manufactured in Denmark by:

