

Mini Themage AYJ-M08

迷你热玛吉 AYJ-M08
Mini Themage



Wrinkles removal
Acne scars
Textural irregularities
Superficial skin lesion

Power Input	AC100-240V
Pweper	AC:9V 500mA
Power	5W
Size	20X4.5X4.5cm
N.W	150g

Thanks for your purchase for our beauty device! Before Using this device, we strongly recommend reading through user manual. Please keep this manual in good condition so that you can look through any time .

Use Manual

Catalogue

Beauty Tips.....	01
Safety Caution	02
Product Structure.....	03
Features	04
Product Effect.....	05
Method Application.....	06
Rejuvenation skill	07
Skin Care.....	08
Beauty Head Replacement	09
Parameter Specification	10
Maintenance.....	11

一: Beauty Tips

A: What is rf ?

R: Radio F:Frequency

RF beauty treatment is a heat penetration through the radio communication apparatus changes in modern FM transmission of electromagnetic energy to the subcutaneous tissue, the dermis to stimulate collagen contraction, regeneration, dermal thickness and density, pulling fill wrinkles, remove scars and restore skin elasticity and gloss; commonly known - Mesotherapy, electric lift.

B: The 3 feature of rf

- 1) Safe non-surgical, non-invasive, no anesthesia, no injection, no pulling pain, noninvasive trace, no recovery high security.
- 2)Quick immediate, tightening of the skin care can make a relaxation of the skin instantly felt a raise, tight treasure effect.
- 3)Hold time long lasting, efficient promote skin collagen continue, 2 to 6 months shallow wrinkles to disappear, restore skin rejuvenation state.

C: Function


Tighten the skin, shrink pores.

Deep Wrinkle facial contours modified to remove a double chin, face-lift.

Activate collagen cells and increase skin elasticity and luster whitening

Fade dark circles eye.

Dilute head pattern.

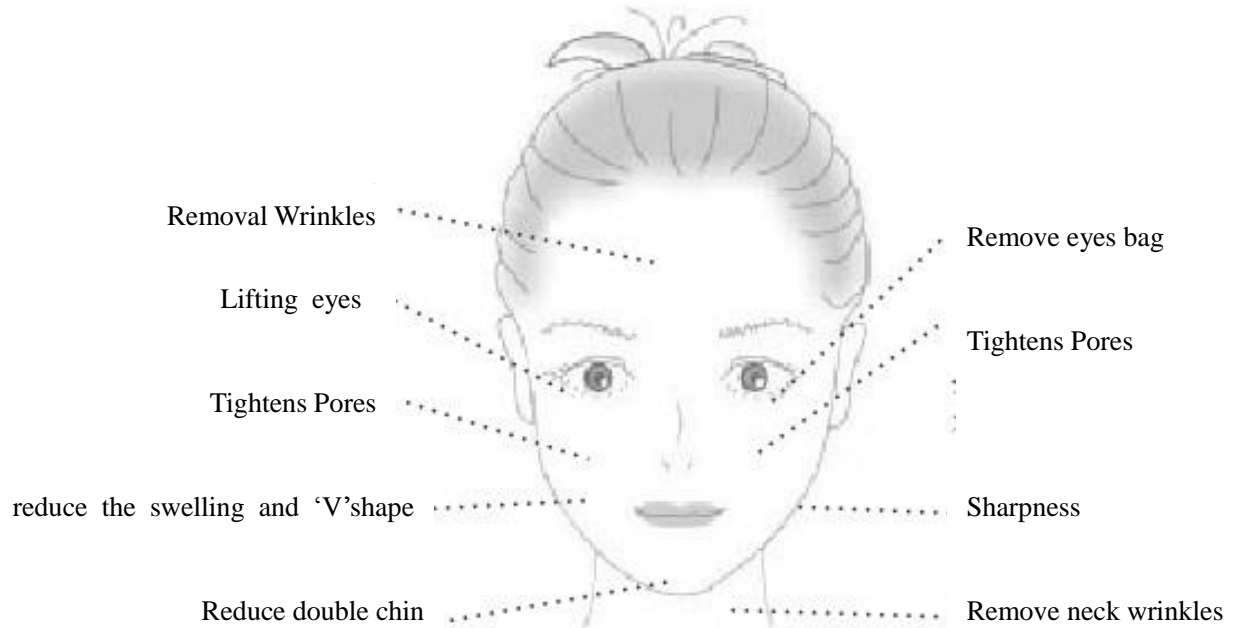

二: Safety Caution

- ★ The instrument righteousness special beauty for healthy skin care equipment, please do not used for purposes other than cosmetic.
- ★ Please observe the operating instructions and use frequency content, initial use or sensitive skin, try to shorten the usage time.
- ★ When the body appears problem , please stop using it.


Attention of the before rf

- 1, if the therapeutic range there have injected hyaluronic acid, collagen within six months , is not recommended.
- 2, if the therapeutic range was injected artificial fillers (such as silica gel), please avoid the area.
- 3, have heart disease, fitted with a pacemaker and pregnant women do not use.
- 4, treatment prohibited him wear any metal objects.
- 5, Do not make-up, when treatment or pre-treatment cleansing facial.

Attention of after RF radio frequency treatment

- 1, please enhance moisture and sun protection.
- 2, three days do not wash with hot water (water temperature can not exceed).
- 3, within three days, please temporarily not hot springs and sauna.

三：Product Structure

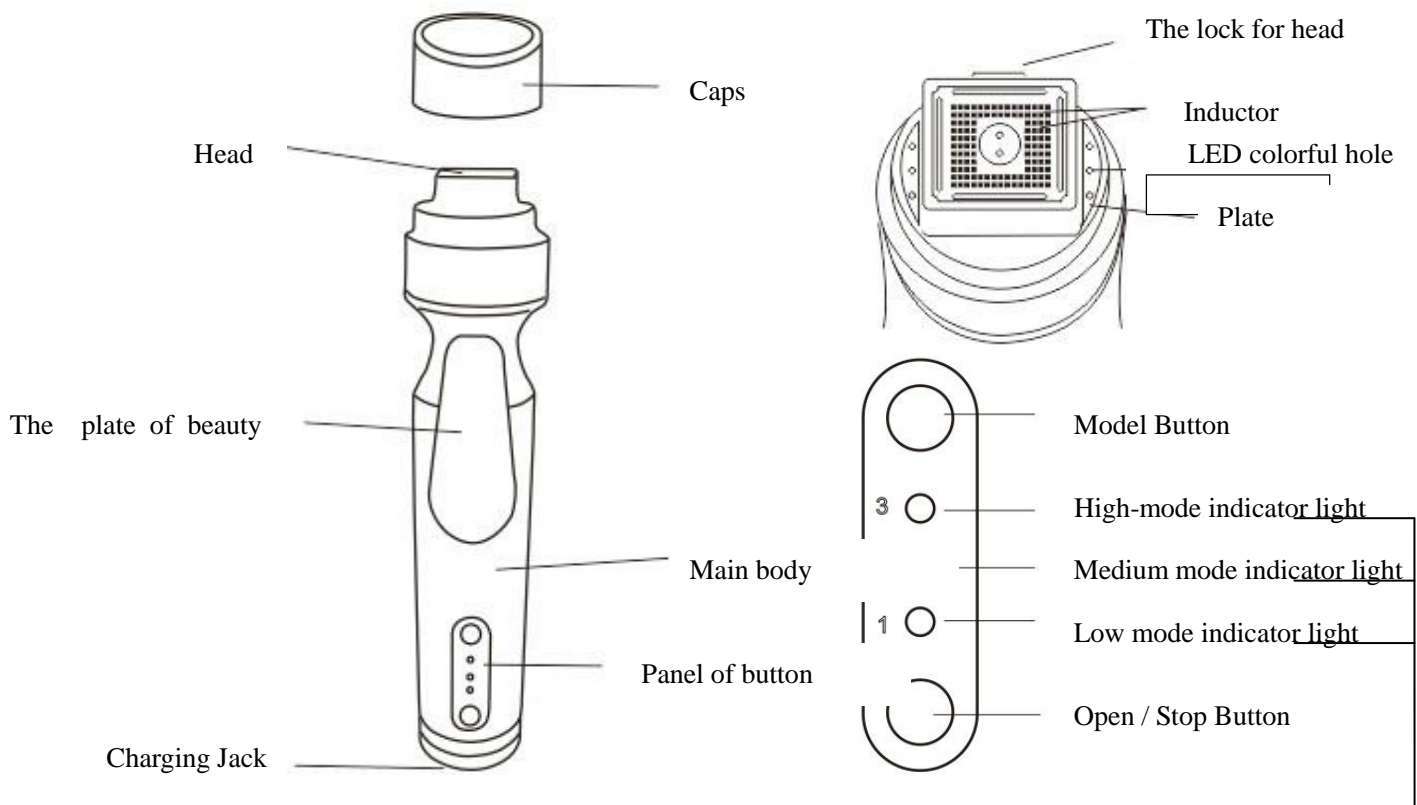

When low voltage is flashes, suggesting the need charging

四：Features

- ◆ dot 144 dot matrix electrodes beauty head, efficient deep repair skin, treating the symptoms, solve your problems and shortcomings!
- ◆ dot electrode beauty head arbitrary replacement design to avoid cross-use, better sanitation and health.
- ◆ once saturated charge (lithium battery), continuous use of three one four hours, 360 unlimited use of space, get rid of the "pull worried about hanging."
- ◆ appearance mini lightweight, portable, anytime, anywhere to meet your care needs.
- ◆ Smart John the detection (cosmetic skin contact output before header automatically open), more secure and energy efficient.
- ◆ Additions infrared phototherapy function, improve blood circulation, promote

metabolism, eliminate facial fatigue.

- ◆ high, medium and low intensity output modes, for various groups and skin types.
- ◆ intelligent low voltage alarm, rechargeable, saturated tips makes the human use function.

五: Product effect

RF Repair skin , tissue

Radio frequency (abbreviated RF) Beauty is a non-hand derivative, a new cosmetic procedure. It uses capacitive coupling electrode to deliver radio frequency energy and produce an electric field across the surface of the skin into the subcutaneous tissue. In the high frequency of 1Mhz and dot electrodes, this field changes per second 86.4 million times or so, in order to restore the electrode rapidly changing, the charged particles in the skin at the same frequency conversion on the direction of the subcutaneous tissue of the natural resistance of movement at this time produce natural energy, dermis collagen in Celsius 45 a 60 degree temperature, it will have an immediate contraction of the features to slack skin in the treatment of complex immediately felt pulling up, tight face lift cracked fruit. When the endless production of collagen, it will make the skin dermis layer thickness and density increase, fill wrinkles, eliminate scars, restore skin dread and gloss, make skin white and smooth.

Thermal effects of radiofrequency waves to increase blood circulation, accelerate metabolism, improve local tissue nutrition, enhance enzyme activity.

Deep import

RF has a strong deep import function; fit contain "active gene repair function" fine completed or skin care products, make deep into care better. In the two together, in order to stimulate the regeneration of skin collagen deep soft tissue reconstruction, newborn skin to achieve whitening, skin, anti the overall effect.




The care of Infrared light

Infrared is an electromagnetic wave of life, the wavelength range of 4 to 14 microns, with strong permeability, can penetrate human skin 15 cm, resonate with human cells and accelerate the generation of new Chen pavilion. When the infrared absorption by the body tissue, it will produce a kind of soothing and beneficial to human body heat, can penetrate into the human group, and repair body cells. The physical properties of infrared, it has: Activate activity of biological macromolecules; promote and improve blood circulation; increase metabolism;

inflammation, swelling and other health care effects.

六: method of application

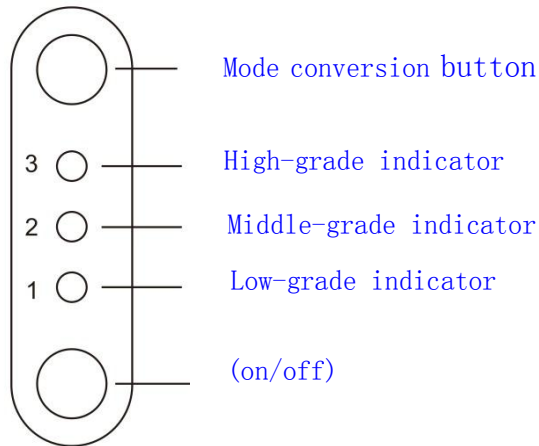

1. Cleaning the skin

2. Take a suitable amount of gel or essence;

apply face, neck and other body parts of nursing

3. Press the button (ON/OFF), start the machine.

4. Press mode selection button (MODE),

adjust the suitable energy level.

(Please try in the arms feeling before operation for the first time, to choose your suitable mode according to what you can accept the degree of stimulation and heat.)

1) Press one time the low-grade indicator flashes

(the green light is on)

2) Press two times the middle-grade indicator flashes.

(the green light is on)

3) Press three times the high-grade indicator flashes.

(the green light is on)

5. According to different nursing position, use the corresponding nursing methods.


6. Infrared light care automatic induction open with nursing technique simultaneously.

7. Press the (ON/OFF) button to close the machine after finishing the nursing.

Skin will feel heat or a slight tingling in the process of nursing. Some users will feel a slight redness of the skin after using. This phenomenon usually disappears within a few minutes after using. If the user feel abnormal skin fever in using, suggest to stop using and let the skin cool.

七: Rejuvenation skill

Put the beauty head close to restored area and stop for a minute.

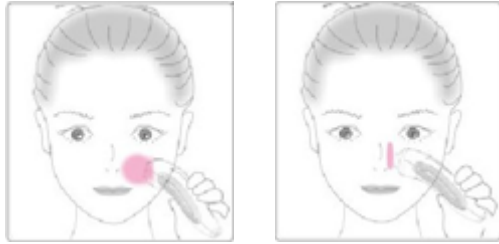

(Suggestion: low-grade mode about 10 seconds, middle-grade mode about 5seconds, high-grade mode about 1 second. Long retention time can cause burns to skin.)
Then move to the other restored area to repeat this action.
This process just as a "stamp" on skin, until finishing the whole nursing.


Special note: avoid eye area, Delicate skin around the eye care please choose low mode.

八: Skin Care (recommend)

A: Dredge lymph


Dredge and awaken lymph, enhance the detoxification and remove edema.

Put the beauty head closely to the collar bone near artery location; make an effort to outgoing circle to behind the ears, and then effectively down along the aorta to the starting point(the collar bone near artery location), repeat the action several times.

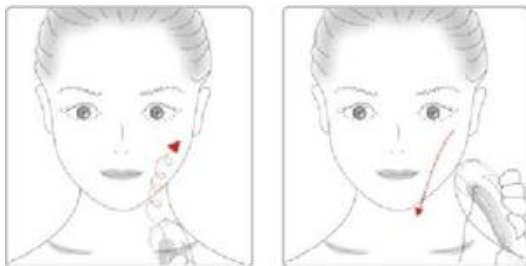

B: Removal double chin, reduce wrinkles

Put the beauty head closely to the double chin, force back to the behind ears, from behind the ears down to the collar bone near the big dolphin position, and then on the contrary, by the clavicle pull up to the jawbone lightly, again forward to double chin position. Repeat the action several times.

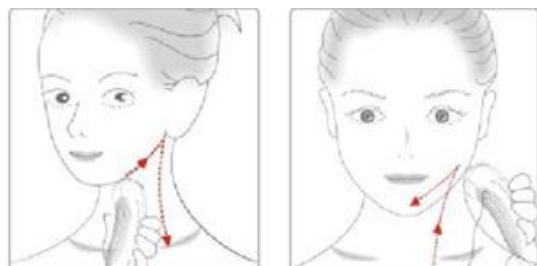


C: Shaping V outline

Put the beauty head closely to chin, circle to earlobe, along the ear again back to the behind ear, from behind the ear down to the collarbone. Then in the opposite direction, from the clavicle pull up to the jawbone lightly, then again forward to chin position, repeat the action several times.

Put the beauty head closely to the mouth, Looped outward to the ear, along the ear again back to the behind ear, from behind the ear down to the collarbone. Then in the opposite direction, from the clavicle pull up to the jawbone lightly, then again forward to mouths position, repeat the action several times.

Put the beauty head closely to the cheekbone, Looped outward to the ear, along the ear again back to the behind ear, from behind the ear down to the collarbone. Then in the opposite direction, from the clavicle pull up to the jawbone lightly, then again forward to cheekbone position, repeat the action several times.

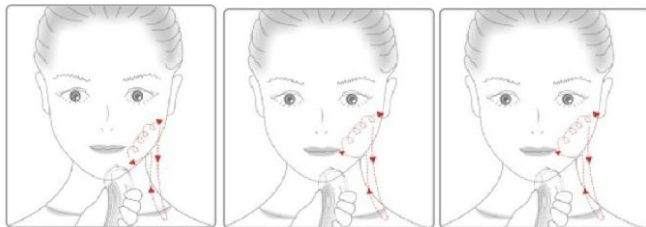

D: Remove wrinkles back to youth

Promote the collagen regeneration, increased flexibility, eliminate wrinkles
Put the beauty head closely to chin, draw the C word hard to philtrum-philtra, then from the philtrum-philtra to the chin, draw the C again lightly. Repeat the action several times


Put the beauty head closely to nasolabial folds, looped outward to the top of the nasolabial folds (beside nose), light force down again to the nasolabial folds. Repeat the action several times.



E: straightening nose and shrinking pore


Promote the collagen regeneration, increased flexibility, shrinking pore
Put the beauty head closely to the nose, outward force on

the ring to the bridge of the nose. Repeat the action several times.

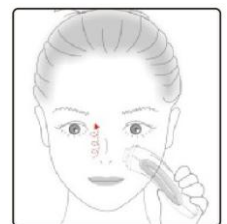


F: Remove eye bags, wrinkle around eye and dark circle

Put the beauty head closely to the eye bag, about 1 cm away from the eyes, pull to temple. Repeat the action several times

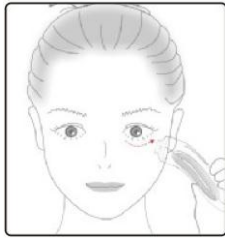

九: Beauty Head Replacement

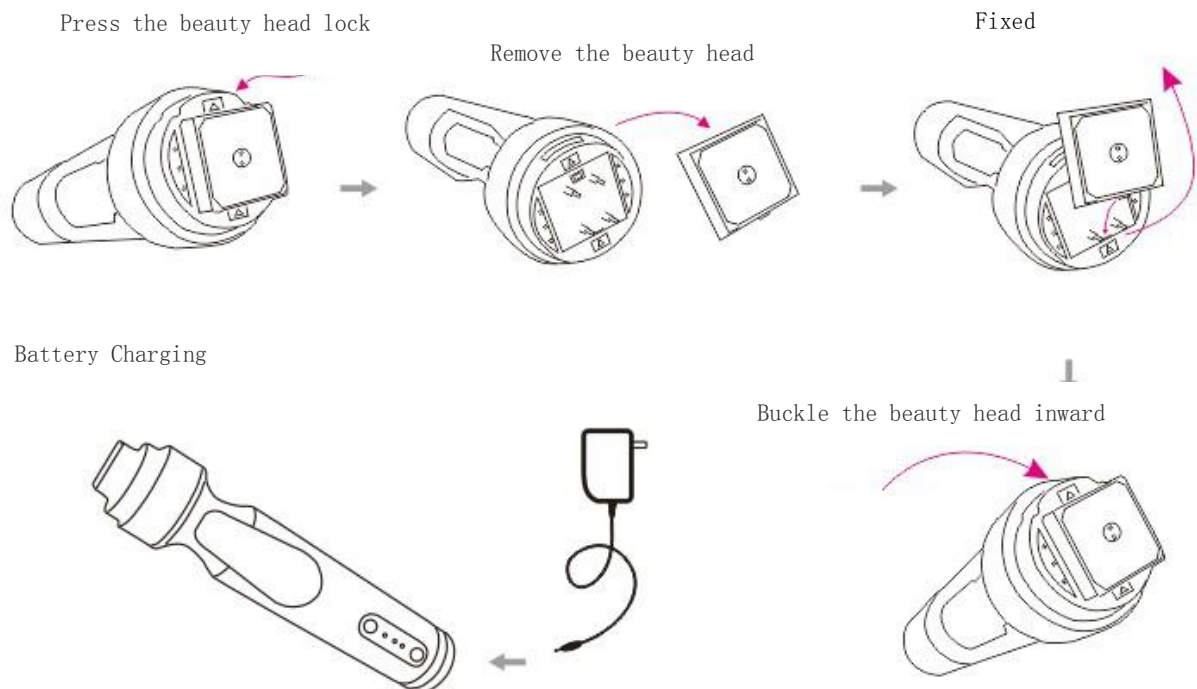

十: Parameter Specification

power input	AC100-240V
power output	AC:9V 500mA
power	5W
size	157×40×42mm
weight	150g

十一: Maintenance

- ◆ Please use a sponge or wet towel to wipe airframe after using; do not wash with water or soaking in water.
- ◆ To clean the beauty head, use sponge to take special pure water or alcohol (70% propylene), wipe the beauty head parallel. If bust beauty gel is used, beauty head should be removed to clean up the whole frame.
- ◆ Do not use sharp objects or metal to clean beauty head, so as to avoid injury of electric circuit.
- ◆ After using place the machine in a cool, dry environment
- ◆ Please charge in time, so that the battery in the fully charged state

Repair Record	
Name :	Phone :
Address :	
Item no :	Number of product :
Purchase date:	Dealers (chapter):