



INFINITE THERAPEUTICS 68 Route 125 Kingston, New Hampshire 03848 U.S.A.

ph: (603) 347-3006 Fx: (603) 642-9291



Thank you for your purchase of this product. Before using, please read this USER MANUAL carefully, and pay special attention to the safety precaution to ensure proper operation and use.

Please keepthis USER MANUAL properly for further reference.



# **CONTENTS**

SAFETY PRECAUTIONS	1
NAME AND FUNCTION OF EACH COMPONENT	3
Exterior structure	3
Interior structure	3
Power supplybox	4
Instruction on each component of main controller	5
Function description on buttons	5
Description on VFD displayer	8
Functions and features	10
PREPARATION BEFORE USING	11
Installation position	11
Installation method	12
Move method	18
Turn on power supply	18
Chair adjustment	18
Before sittingdown	18
User selection	18
PRODUCT STORAGE AND MAINTENANCE	19
TROUBLESHOOTING	20
TECHNICAL SPECIFICATIONS	21

# TECHNICAL SPECIFICATIONS

Product name: Infinity Massage Chair

Product Model: 8000 Series
Rated voltage: 220-240V~
Rated frequency: 50/60 Hz

Rated power: 260W Safety design: Class I

(21)—



# **TROUBLESHOOTING**

No.	Problem	Possible Cause	Solution
1	You hear a whiring or rubbing sound.	It is anormal sound made by the working motor or the sound is produced when the massage ball is rubbed against the cloth cover.	It is normaland no need fortreatment.
		The powersupply plug is detached.	Insert the power supply plug into the socket again.
2	The machine can not be started up.	The powersupply switchis not turned on.	Turn onthe switch
		The fuse of massage chair blew out.	Replace a new one with the same mode
3	Some music can not be played normally	To checkwhether the format of musicis MP3.	Replace the MP3 music file
4 heard w	Alarm sound is	The U-disk is not connected well.	Insert the U-disk again.
	heard when pressing the music	The U-disk isinitializing	Wait 5-10 seconds
	control button.	The format of music files in U-disk is not unmatched.	To check MP3 music files with the help of computer.
No sound f the speake earphone	No sound from	The volume is turned to lowest.	Turn thevolume up
		The waiting time of the massage chair is over 2 hours.	Press the <b>(b)</b> button again
6	Backrest or legrest can not ascend or descend.	Over pressureis put on backrest orlegrest.	Reduce the pressure on backrest or legrest.
7	The airbag installed at seat cushion or legrest can not work.	To check whether there is air tube blockage.	To clear the airtube
8	Some fault code ( " E0-E7") is displayed while the main controller is buzzing.	The massage chair has automatically detected some exception.	Cut off the powersupply and start up themachine again to check whetherit can return to normal, if not, please send it to the professional maintenance center for repairand maintenance.

If abnormal phenomenon still remains afterabove measures have been taken to the massage chair, please contact local distributor for repair and

Note: Do not dismantle the product by yourself, and our company shall shoulder no responsibility to the damage of the product or the injury of people caused by improperuse! Thank you for your cooperation.

# SAFETY PRECAUTIONS

## Safety Precautions

- ▶ Before using, please readthis "SAFETY PRECAUTIONS" carefully to ensure proper operation.
- ▶ The following "SAFETY PRECAUTIONS" is of great help to your correctly use of this massage chair.
- ► To clearly indicate the degrees of danger and damage, the classes of accidents caused by incorrectuse are classified into "WARNING" and "CAUTION", instructions marked as followings should be strictly followed:

<b>MARNING</b>	Improper operation may cause injury to users.
<b>A</b> CAUTION	Improper operation may cause users in jury or product damage.

Forbidden. Specific instruction will be indicated in the diagram drawings. (e.g, the left drawings means disassembly is not allowed.)
 Compulsory Specific instruction will be indicated in the diagram drawings. (e.g., the left drawings means to disconnect the power plug)



	We recommend that those with any of the following conditions should consult a doctor before using.			
COMPULSARY	Those whoare using electronic medical instrument imbeddedin the body, such ascardiac pacemaker etc.     Those thatare being treated by doctors, especially those who feel uncomfortable.     Those whoare malignant tumor patients, heart disease patients, acute patients.	<ul> <li>4. Womenwho are pregnant or menstruating.</li> <li>5. Those withosteoporosis or fracture spine.</li> <li>6. Those withskin disease or those whose skinsare injured.</li> <li>7. Those whosebody temperature isover 100 degrees F. (Febrile Phase).</li> </ul>		
	Do not allow the mentally disabled to use the chair alone. (It is OK under supervision)			
	Do notallow children play on massage chair.			
	Do not allow more than one person to use this massage chair at a time.			
PROHIBITED	If the powercord or power plug becomes damaged in any way you must disconnect from the main supply immediately and consult an authorized service engineer for repair or replacement.			
	Do not damage, process, excessively bend, pull hard, twist or knot power cord.			
	Juveniles are not allowed to use this product.			
NO DISASSEMBLY	Anyone else, except maintenance personnel specified by our company, should not disassemble, maintain, or remodel this massage chair. If any problem, please contact local dealer immediately.			

## SAFETY PRECAUTIONS



## CAUTION

Please do not continuously use the massage chair over 20 minutes, or limit a massage to a particular part of body for more than 5 minutes.

Please put the massage chairon flat ground for use.

When adjusting the legrest or backrest, please do not exert excessive force on it; in addition, keep children and pets well clear.

When removing power plug, please hold thein sulation plug, do not pull power cord.

Please do not use the power supply that is beyond the specified electric voltage, as this may cause fire disaster or damage massage chair.



Before using, please uncover the back cushion to check whether backrest or fabric is damaged; if damaged, please stop using this massage chair immediately, and contact franchiser for repair immediately.

Please do not operate the controller with damp hands, and do not put heavy objects on controller.

Weight limit is 264 lbs.

Please do not use other therapeutic equipments together with this massage chair.

During use, if you feel uncomfortable, please stopuse immediately and consultyour doctor.

This massage chair is designed and intended for non commercial household use only, do not use the massage chair for medical treatment.

The massage chair is equipped with heating function, so people who are insensitive to heat must use it carefully.

Always switch off power and remove the plug from themains socket immediately after useand before cleaning.

To avoid causing discomfort, do not use the massage chair immediately after meals.

Do not fall into sleep on massage chair, when the massage chair is working.



Do not use the massage chair when you are drunk or feel bad.

Do not sit on backrest, legrest, or armrest to avoid accidental injury or machine fault.

If the upholstery becomes dirty, please use onlyknown good quality leather cleaner, do not use chemical cleaners, petrolor detergents

Please do not smoke while using the massage chair.

etc, and do not spray pesticide onit.



Before maintenance, please pull outpower plug. Inaddition, please do not insert or unplug the power plug with damp hands, to avoid electric shock.

**REMOVE PLUG** 

When the massage chair is not in use, pleases witch off power and remove the plugfrom socket. When power is off, please cutoff power supply immediately and pull out the powerplug from socket.

## PRODUCT STORAGE AND MAINTENANCE

#### Maintenance

- ▶ Please twist the cotton cloth which has been dipped in neutral detergent before wiping the backcover, hand rest, legrest, air tube and some other plastic components, then use a dry cotton cloth to wipe them again.
- ▶ Use dry cloth to wipe the controller, DVD, and power supply box.
- ▶ Use a piece of cloth that has been dipped with neutral detergent to wipe the back cushion and pillow, then use a piece of wet cloth to wipe them, and finally putthem to a well-ventilated place for natural drying.



### Caution

- ▶ Make sure to pull out the power supply plug before maintenance. Do not use wet hands to pull out the plug to avoid electric shock. When cleaning the back cover, please be cautious of detergent enteringinto the power supply box!
- ▶ Do not use abrasive cleaning agents like benzene or thinner to clean the product, otherwise, the color of upholstery may fade or the massage chair may crack or be eroded.
- ▶ Do notiron the upholstery.

## Storage method

▶ Wipe off dust and dirt from the product after use, if the product is not to be used for a long period of time, please cover the product with a dust cover (selfprovided).



## **A** Caution

▶ Do not store the product under direct sunlightor in places of high temperature. This may cause the color of the upholstery to fade.

### Gentle reminder

During daily use, please check whether the following symptoms exist:

- ▶ Scorched smell exists.
- ▶ When handling the power cord, the powerdoes not cycle on then off.
- ▶ The power supply cord is hot.
- ► Any otherther abnormal symptoms




### Caution

- ▶ Please stopusing the massage chair to avoid breakdowns and accidents from occurring if any of the above symptoms exists.
- Immediately turn off the powersupply switch and pull out the power supply plug, then contact the local distributors or headquarters for repair and maintenance.
- ▶ Except for normal operation and maintenance as mentioned above, all other repairs and maintenance must be diagnosed and performed by our authorized dealer or service agent.

## PREPARATION BEFORE USING

## Machine Adjustment

- ▶ Put the backrest to a limiting position, raise the legrest to the highest position and stretch it out to the longest position with the controller, meanwhile, pleasemake sure sufficient space is left for the massage chair.
- ▶ After testing each function, please reset the chair, and turn off the power supply switch. (Please refer to Buttons function on Page 5 and VFD display on Page 8 while you are choosing massage program)

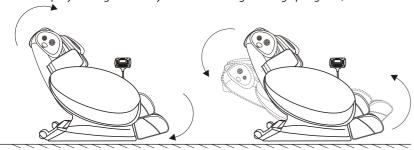

## Before sitting down

- ▶ Make sure the leg massager is well connected with the main chair.
- ▶ Before sitting down, make sure all the components are back to their original position, if not, pleaseturn on the power supply switch and reset the components.



## Caution

- ▶ Do not sit on the massage chair when the legrest is raising.
- ▶ Do not put heavy burden onto the legrest suddenly, otherwise, the mechanical structure may be damaged.


## User selection

Based on your physical condition and the massage time you require, you can choose any of the following automatic procedures to enjoy targeted service.


- ▶ Recovery: It is suitable for doing it after work or one hour later after supper to help people recover their spirit.
- ► Extend: It is suitable for officeworkers or drivers.
- ▶ Relax: It is a kind of leisure massage; especially suit for morning or noon relaxation.
- ▶ Refresh: It is suitable for alleviating the sore or pain at back, waist or neck. It can massage at a particular sore or pain part of body with semiautomatic functions.

## NAME AND FUNCTION OF EACH COMPONENT


### Exterior structure

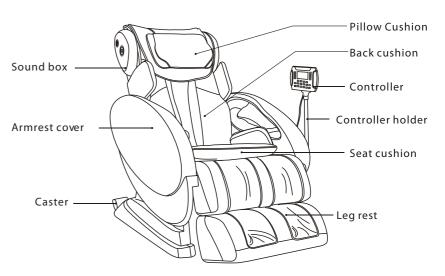

### Interior structure

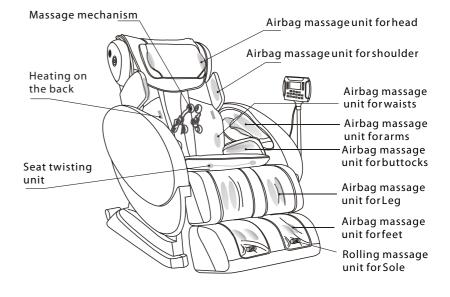

# Power supply box

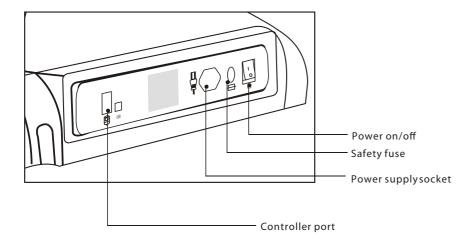




## A Caution:

- ▶ Do not let children play on the chair in of accident happens due to the power supply box.
- ▶ Please power off the chair and unplug the chair from the power supply socket after the massage finishes
- ▶ Do not use the chair during thunderstorms.
- ▶ Please power off the chair and unplug the chair from the power supply socket when the chair is not in use, so that children won't accidentally get the chair started.
- ▶ Please storethe chairif it is not in use for a long time.
- ▶ Please do not move the massage chair on easily damaged surfaces, such as wooden floor and carpet.
- ▶ Do not use the massage chair in damp environment to avoid electric leakage or electric shock.
- ▶ Do not place the massage chair near stove or where it will be exposed to strong sunlightfor long periods as this may damage the leather and bleach the color.

## PREPARATION BEFORE USING

#### Move method

- ▶ The equipment can be easily moved with the caster installed at the bottom of the product.
- ► Ascend the backrest to the highest position.
- ► Turn off the powersupply switch and pull out the plug.
- ▶ Raise the front part of the chair as shown on the picture below, and move the chair to a desirable place.

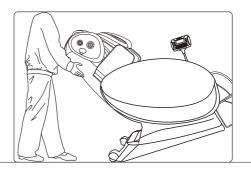




# **A** Caution

- ▶ Do not move the chair while people are sitting on the chair.
- ▶ When moving the massage chair, please hold its base and backrest instead of its armrests (at least 4 people).

# Turn on power supply

- ▶ Put the power supply plug into a well grounding three-corepower supply socket.
- ► Turn on the power supply switch located at the bottom of backrest.

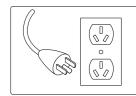


# Caution

- ▶ Before starting up the machine, please make sure the power supply cordand exposed electric wire have not been damaged.
- ▶ Before putting the plug into the power supply socket, please make surethe powersupply switch is turned off.

## PREPARATION BEFORE USING


### Step 6 Installation of back cushion and buffering cushion

► Fix the buffering cushion and backcushion on the backrest with zipper, and install theheater accessory on the backcushion, connect the connector of pillow and waist's air tube as the illustration shows. Connect the lower part of the cushion with seat cushion by linking tape.

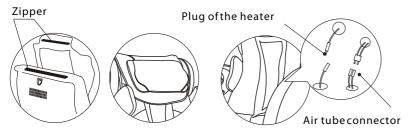


### Step 7 Installation of seat cushion

- ▶ Place the seat cushion on the seat shelf, and fix the air tube.
- ► Connect clothon the seat cushion with the linking tape.

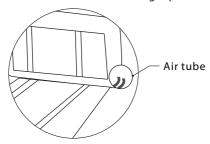

#### Step 8 Other accessories

- Connect the controller, wire to respective jackon the power supply box.
- ▶ Place the controller in the tray.
- Connect the power and switch on the chair, use the controller to return the back rest to normal position.

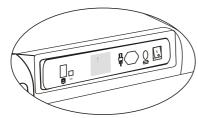

## NAME AND FUNCTION OF EACH COMPONENT

## Instruction on each component of main controller



## Function description on buttons



When the massage program is in POINT mode, keep Upor Down button pressed to pinpoint massage position.



▶ When the body sensing is on and the indicating lightfor shoulder adjustment is on, press Up and Down button to pinpoint shoulder position.



- ▶ Press this button to turn on or off the waist twist massage, and the relevanticon will brighten or go out.
- ► When waisttwist massage is underworking, press this key to choose proper intensity within 1-3 degree. Under maxintensity degree, press this key to stop this function.

  (Notice: When the icon is flashing, it is the fault prompt.)



Press thiskey toturn on/off heating function on backand waist.



► Keep pressing this key to adjust backrest and legrest to a desirable position.



► Keep pressingthis keyto adjust legrest to a desirable position.

## NAME AND FUNCTION OF EACH COMPONENT











Pause/Play

Stop

Previous

Next Volume-

Volume+

Note: press Previous once to return to the beginning of the present episode, continuously press this button twice to return to previous episode.

## ZERO GRAVITY



▶ When this icon is constant on, Zero Gravity starts. Zero Gravity includes weakand strong status, which cater for different people needs. When you feel the massage force is weak and need stronger massage, please press Zero Gravity button again for strong status. When it reaches strong zero gravity, press this button again to cancel zero gravity. The chair returns to original position when it is in regular massage mode and the Zero Gravityicon is off.



#### Recover: Recover autofunction on

- ▶ When this icon is off, pressthis button to start thewhole auto program, which is mainly forstrong massage on big scope.
- ▶ When the icon is constant on, the chair is in Recover mode.



#### Refresh: Refresh autofunction on

- ▶ When this icon is off, press this button to start the whole auto program, which is mainly for moderate strength massage on big scope.
- ▶ When the icon is constant on, the chair is in Refresh mode. (Remark: This program has stretch massage.)



### Power button: Massage chairon/off

- ▶ When this icon is flashing, all massage functions are not available. Press this button to start massage, and then all massage function buttons are flashing for choice. Press anyone of these to start massage, whose defaulted time is 20mins.
- ▶ When this icon is constant on, press this button to end all massage functions, the chair goes back to its original position.
- ▶ When this icon is off, all buttons on the controller are not available.



#### Extend: Extend auto function on

- ▶ When this icon is off, pressthis button to start thewhole auto program, which is mainly forweak massage on big scope.
- ▶ When the icon is constant on, the chair is in Extend mode.



#### Relax: Relax auto function on

- ▶ When this icon is off, press this button to start the whole auto program, which is mainly for moderate strength massage on whole body.
- ▶ When theicon is constant on, the chair is in Relax mode. (Remark: This program has air pressure massage of pillow.)

#### **AIR BAGS**




- ▶ When this iconis off, and massage function is already on, press this button to start airbags massage mode.
- ▶ When this iconis constant on press this buttonto end airbags
  - (Remark: In automatic program, the airbagsand roller is linked.)

## PREPARATION BEFORE USING


#### Step 4Installation of the speaker

▶ Open the screw cap cover on the speaker; insert the line from upper side of backrest to the line from the speaker, then aim the speaker to the hole on upper side of the backrest, and fasten them by screws. (It is the same installation method for the other side)

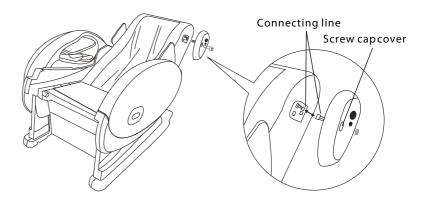

#### Step 5 Installation of controller holder


- ► Take off the small cover on the left armrest
- ▶ Aim the installation hole of the controller holder to the hole on armrest, and thenfasten themtogether by screws, lastly place the small cover.

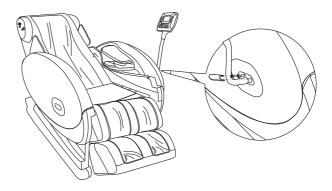


## PREPARATION BEFORE USING

Step 3 Installation of learnassage shelf

- ▶ Insert the leg connection wire plug to its jack.
- ▶ Connect air tube from the leg massage shelf with its port on the seat
- ▶ Aim the coupling shaft port on the leg massage shelf to the coupling shaft on the seat and fasten them with one screw each side.
- ▶ Paste the folding PU cover on the leg massage shelf with linking tape on the cloth of the seat.

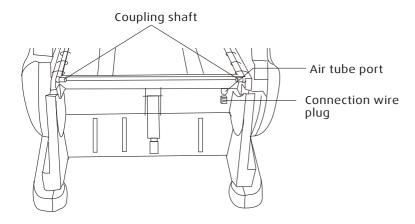

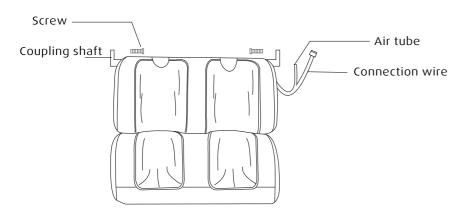



## NAME AND FUNCTION OF EACH COMPONENT

**SOLES** Press this key to turn on/off rolling massage function on sole



▶ When the indicator light is off, press this key, the light will turn to flash status, then interval rolling massage on sole is started and press twice, the indicator light will be constantly on, then continuous rolling massage on sole is started, press it the third time, rolling massage on sole will be turned off. When the indicator is flashing or constantly on, the Product is under rolling massage on sole (Remark: In automatic program, the airbag and roller is linked.)



▶ When the massage functions are under manual mode, the indicator light is constantly on, this key is operable, press it to circularly choose full, partial and point, and the back massage mechanism correspondingly conductfull, partial and point mode massage.

▶ When the indicator light is off, this key is not operable



▶ When the massage functions are under manual mode, the indicator light is constantly on, this key is operable, press it to circularly choose narrow, medium and wide, and width of the crank arm correspondingly change.

▶ When the indicator light is off, this key is not operable.




▶ When the massage functions are under manual mode, the indicator light is constantly on, this key is operable, press it to circularly choose 6 different speeds.

▶ When the indicator light is off, this key is not operable.



▶ When the indicator light is off, if massage function is under working.

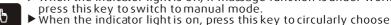

kneading, pummel, kneading \$knocking, knocking, shiatsu and rhythm massage.







Arms & shoulders, upper body, legs & feet

▶ When massage function is under manual mode, if the indicator light is off, press this key to open turn on the corresponding function, then the indicator liaht will be on. ► When the indicator light is constantly on, press this key to stop corresponding airbag massage function.



▶ When the indicator light is constantly on, press this key to choose strength among weak, medium and strong.

▶ When the indicator light is off, this key is not operable.



# **Caution:**


▶ Every auto program has lots of massage action on shoulder. The chair conducts body shape sensing on the initiation of auto program (or after electric cylinder position is adjusted), so the userhas to sit wellon the chair at this moment until body shape sensing finishes. (i.e. shoulder adjustment on VFD display stops flashing).

▶ To better adapt to the environment, the controller of this chair has dynamic sampling. Within 5 seconds after the chair powers on, please don't use hand or other metallic substance to touch controller display; otherwise it will affect its button sensitivity.



## NAME AND FUNCTION OF EACH COMPONENT

## VFD display description

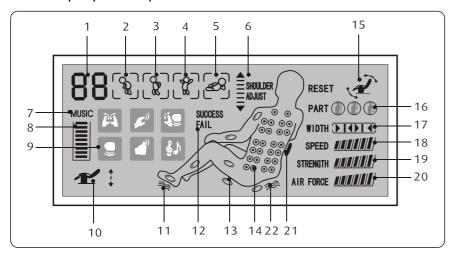

- 1.Time display: It shows working time upon the start of massage with defaulted time of 20 minutes.
- 2.Recovery indicator
- 3.Refresh indicator.
- 4.Extend indicator
- 5.Relax indicator

Under the massage modes: Recovery, Refresh, Extend or Relax, the icon will be on, but not under other modes.


- 6.Shoulder position adjustment indicator
  - During body shape sensing, after it automatically senses on the back part of shoulder, this area is lighted to remind the user to micro adjust shoulder poison by Shoulder Up/down adjustment
- 7.USB Flash disk connection indicator
- 8.MP3 Volume indicator
- 9. Massage modes: to show present massage mode

Kneading u ,tapping u ,sync u ,knocking u ,shiatsu u , music u


- 10.Zero Gravity cylinder performance indicator: lie upwards; return downwards.
- 11. Rolling masage on sole indication
- 12. Body sensing status indication
- 13. Part of airbag massage dynamical indiction
- 14. Massage dynamical indication

## PREPARATION BEFORE USING

▶ Make sure the 4 connection sticks completely inserted into the eccentric joints, then use cross screwdriver to fasten the armrest by turning 180 ° clockwise.( The sharp point on the eccentric joint shall point downward " ) It is the same installation method for the other armrest.

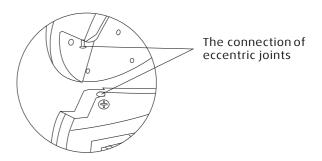

▶ The connection of armrest locks with backrest. Hold up the rear part of the armrest, push inward to let the gap of armrest lock aim at armrest spindle, then push the armrest to fasten them together.(Press downward the handle of the lock to unlock armrest and backrest)

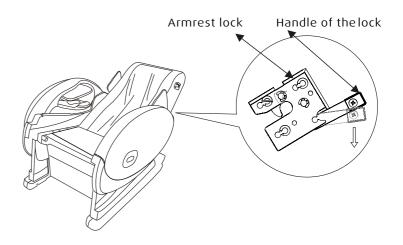




## Installation method

Step 1: Take out the components of massage chair

- ▶ Open the main carton, and take out the quilting product packages, accessory carton and the main body of the massage chair carefully.
- Check whether the components are complete (please check one by one according to the packing list)


## Step 2: Install left and right armrests

▶ connect air tube port under the armrest to its corresponding output air tube. (It is the same way for the other side)

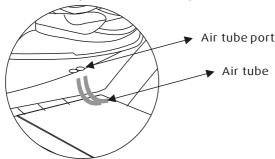

▶ The connection of armrest and armrest shelf: Aim the three connection sticks of armrest to the three eccentric joints on the armrest shelf, then press downwards on the armrest and fasten it by the eccentric joint.

(If the rear screws are obscured by the back cover, and can't be accessed, power on the chair to raise the back.) It is the same installation for the other armrest.

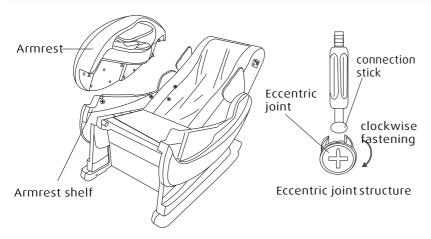





## NAME AND FUNCTION OF EACH COMPONENT

- 15. working indicator of reclining motor for backrest and legrest
- 16. back massage position mode:



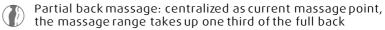

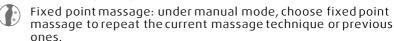

17. back massage width indication3 width adjustable: narrow, medium, wide



18. back massage speed indication 6 degrees adjustable



19. Vibration(twist) strength indication3 degrees adjustable



20. Airbag intensity indication 5 degrees adjustable



- 21. back heating indication
- 22. seat twist indication

# F

## **FUNCTIONS AND FEATURES**

### 1.shouler position

shoulder detecting function can locate exact position of shoulders of different people, thus to ensure the most effective massage on shoulder and neck.

#### 2.Zerogravity

Abiding and imitating principle of space capsule, the Product can give fully relax and most comfortable massage posture, it conforms to ergonomics and keep back, but tocks and legs in a weight less status, thus to ensure people on the Product lay in a most relaxed reclining angle, then weight of body will fully be released on the Products, muscles will also be in a completely relax condition, from which people can enjoy weight less massage experience.

### 3. Four-wheel massage mechanism

Innovated designof four wheels massage mechanism and walking structures, it conforms to ergonomics.

#### 4.6 massage strengths

Under manual mode & anto mode, 6 massage strengths are optional, people can choose the most desirable one according to health condition and massage time.

## 5.airbag pressure massage for arms

upper and lower airbags design to realize wavelike massage on arms.

#### 6. Shoulder massage

Side airbags installed on both sides of shoulders, it can give air pressue massage on shoulders.

#### 7.airpressure massagepillow

Pillow airbags can conduct swinging and kneading massage on neck.

#### 8. Waisttwist and airbag massage on seat

Waist twistfunction in seat, airbag pressure massage on thigh and waists.

9.Airbags pressure massage on legand feet, rolling massage onsole Multi-layer airbag design to ensure quick and abundant massage on leg &feet, three rows of rollers located under sole to conduct shiats u massage on sole, archand hell offeet.

#### 10.Mp3 & DVD (optional)

Mp3 function, DVD is optional as per requirement of customer.

#### 11.Touchkey

Controller is adopt to touch key, it looks more fashional.


## PREPARATION BEFORE USING

### Installation position

Safety distance for the periphery of massage chair

- ▶ The horizontalinterval ofbackrest between the highest and lowest
- points should be over 50cm.

The horizontal interval of legrest between the highest and lowest points should be over 40 cm.

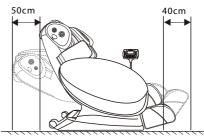

## Floor protection

► The floor maybe worn by the massage chair, so it is advised to put some cushion material (for instance: a carpet) under the massage chair to protect the floor.



# Surrounding environment



Keep themassage chair away from damp places (for example: bathroom) to avoid electric leakage.



 $Keep\ the\ massage\ chair away\ from\ high\ temperature sources\ (for\ example:\ stove) to\ avoid\ fire\ or\ damage\ to\ leather\ material.$ 

## Earthing

▶ This massage chair is a Level One electrical appliance, make sure to use three-core power supply socket and to guarantee the ground wire connected to the power supply socket has obtained good ground connection to avoid electric leakage, electric shock and some other negative effect during use.