




User Manual



Contents

| | |
|---|----|
| Safety Instructions | 4 |
| Meet no!no! HAIR LITE | 5 |
| Before you no!no! | 6 |
| Get ready to no!no! | 7 |
| Working with no!no! HAIR LITE | 8 |
| Thermicon™ Tip | 10 |
| Maintenance | 11 |
| Troubleshooting | 12 |
| Technical Specifications | 13 |

Safety Instructions

WARNINGS – to reduce the risk of burns or injury

- Do not use this device on the genitals, nipples or face.
- Do not use this device on sunburned skin, open cuts, bruises, blisters, moles or on raised areas such as scars.
- If you have any skin condition other than those listed above, like cold sores, herpes, severe acne, please consult your physician before using this device.
- Do not use deodorants, lubricants, liquids, creams, gels, alcohol, water, etc. before use.
- Do not allow the thermodynamic wire to touch your skin.
- Only use this device for its intended purpose as described in this manual.
- Do not leave the device unattended while plugged in.
- Close supervision is necessary for children under 14.

CAUTIONS - to reduce the risk of fire, electric shock, etc.

- Do not reach for a device that has fallen into water while it is plugged in.
- Unplug the power supply from the device before removing it from the water.
- Do not place or store device near a sink or bathtub.
- This device should never be left unattended when plugged in.
- Do not use the Thermicon™ Tip if the wire is bent or damaged in any way.
- Turn off the device while changing the Thermicon™ Tip.
- Do not use the power supply if the cord or plug are damaged in any way.
- Do not use this device near or rest it on flammable objects such as paper or cloth, etc.
- Always unplug the power cord from the electrical outlet after use and before cleaning or maintenance.
- Do not place this device near water or other liquids.
- This device should only be used with the supplied power supply.
- Always turn off this device when you have finished using it.
- Do not open the outer casing of the unit. This might harm the unit and it will void any warranty.

Meet no!no! HAIR LITE

Congratulations, you are about to enter the world of virtually painless hair removal that gives you salon-like results. As you do, it's important to remember that no!no! HAIR LITE is an ongoing hair removal treatment for professional salon-like results at home. This is science, not magic. Results take time, patience and a strong commitment to the treatment's unique procedures. Since every body is unique, results will vary. The average user should notice results in 4-6 weeks, for others it might take a bit longer.

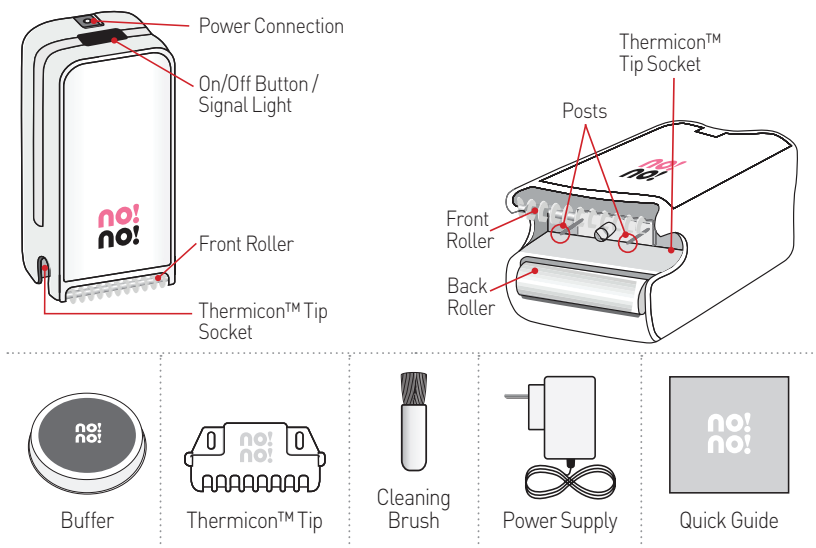


Figure 1

Before You no!no!

Get the best results and performance by learning as much as you can before you no!no! with no!no! HAIR LITE. Below is a brief description about what to expect and what to look out for and how each of the parts work.

- **Test run** - Before you use no!no! HAIR LITE the first time, and before treating new areas, always test no!no! HAIR LITE on a small section.
- **Track your progress!** - Before you use no!no! HAIR LITE the first time, study the area you want to treat. You might even want to take pictures or choose a small area and count the hairs. Then, once every 2 weeks, check your progress. After 4-6 weeks, you should be able to see initial results.
- **The smell that tells you it's working** - no!no! HAIR LITE burns hair and burnt hair smells. But that's how you know it's working. An after treatment cream or moisturizer like no!no! Smooth After Treatment cream will help mask the smell.
- **Getting used to the sensation** - Some people experience a mild tingling sensation when they first use no!no! HAIR LITE, especially around sensitive areas. This sensation will fade away completely after a few treatments.
- **For the sensitive type** - If your skin feels uncomfortable during treatment, take short breaks after repeatedly gliding on any one area. If it continues, stop treatment completely.
- **Buff, the finishing touch after you no!no!** - It is essential that you use the Buffer to buff away the crystallized hair that remains after treatment. Buffing exfoliates your skin and leaves your skin feeling silky smooth.
- **Keep it smooth** - When you finish your treatment, apply a non-alcoholic lotion or cream like no!no! Smooth After Treatment cream to keep your skin moisturized.
- **Patience is a virtue** - Initial results take an average of 4-6 weeks to appear.
- **The first 4-6 weeks** - For best results, use no!no! HAIR LITE at least 2-3 times a week for the first 4-6 weeks. You don't have to wait for a lot of hair to grow back. Just 1 mm of stubble is enough to get the smooth results you are looking for.
- **After the first 4-6 weeks** - You may have noticed that your hair isn't growing in as fast or as thick. You can reduce the treatments to once or twice a week or use as needed to maintain your results.
- **Check the Thermicon™ Tip** - Always make sure the thermodynamic wire is completely intact before use. If it is bent or torn, replace it before your next treatment.
- **Clean the Thermicon™ Tip** - Use the cleaning brush frequently to remove hair debris, especially when working on long hair. This will help prolong the life of your Thermicon™ Tip.
- **Stay Smart!** - Do not use no!no! HAIR LITE on the face, genitals or on and around the nipple area.
- **Stay protected!** Your no!no! HAIR LITE warranty lasts for one year on parts and service. Keep your receipt.

Get Ready To no!no!

1. Your no!no! HAIR LITE comes to you preassembled, with the Thermicon™ Tip already inserted. Simply remove the safety seal from the Thermicon™ Tip socket.

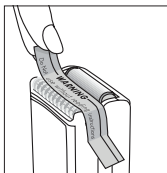


Figure 2: Remove the Seal

2. Always check the Thermicon™ Tip to make sure it is intact. Never use a tip if the wire is bent or torn.

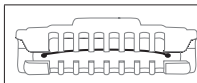


Figure 3:

Bent thermodynamic wire

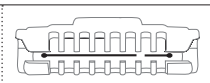


Figure 4:

Torn thermodynamic wire

* For more information, see the Thermicon™ Tip section.

3. Connect the power cord to the no!no! HAIR LITE and then to the wall socket.

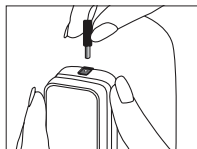


Figure 5: Power Supply Connection

4. Turn the unit on by pushing the on/off button found on the top of the unit, next to the power cord connection. You will need to press for a few seconds until the blue signal light turns on.

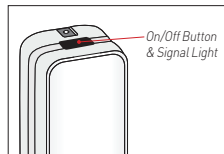


Figure 6: Turn on the button on your no!no! HAIR LITE

*The blue light will flicker until you start gliding.

Working With no!no! HAIR LITE

Skin Preparation

Thoroughly clean and dry your skin. Do not use no!no! HAIR LITE on any area where there might be residual gels, oils, lotions, deodorant, etc.

Mastering the Slow Glide

- Hold no!no! HAIR LITE at a 90° angle with the rollers flat against the skin.

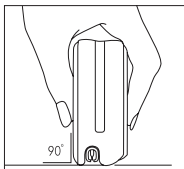


Figure 7

- Gently glide no!no! HAIR LITE over each section and against the grain of hair growth (Figure 8).
- Always glide in a slow, steady, smooth, even motion in one direction. You may need to go over each area more than once, but DO NOT stop in one place.



Figure 8: Glide

- Keep the blue light steady. If it flickers, you may be going too slow, adjust your speed accordingly.
- Do not press down too hard with the device.

Mastering the Buffer

You've tried everything but your skin still doesn't feel smooth? Did you remember to buff? Without buffing, crystallized hair remains in place, giving your skin a prickly feel. Buffing is the essential final step to a no!no! HAIR LITE treatment. Buffing exfoliates your skin, removes crystallized hair and leaves your skin feeling silky smooth.

To Buff:

1. Hold the Buffer by the flaps.
2. Firmly rub the buffer in either a circular or back and forth motion over the treated area.

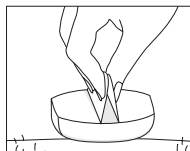


Figure 9: Proper Buffing

Sensitive, Curved and Round Areas

Because of curves and angles, rounded areas like elbows, knees, wrists, underarms and the bikini line require special care and experience. Once you have mastered the Slow Glide on your legs or forearms, you can try no!no! HAIR LITE on more difficult areas.

To work on rounded areas, go slowly and allow no!no! HAIR LITE to follow the natural curve of the body.

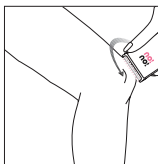


Figure 10: Follow the Natural Curve of the Body

NOTE: The bikini line is more sensitive than other areas. If it becomes uncomfortably warm during treatment, move away from it and work on a different area for a while. Do not go past the bikini line.

Watch the Light

A blue signal light appears from the on/off button as soon as you turn on your no!no! HAIR LITE. It will flicker if you are not gliding at the right speed.

The blue light will turn red when it's time to replace the Thermicon™ Tip. Simply remove the Thermicon™ Tip and replace it with a new one. (See the section on Thermicon™ Tips for more information.)

To Summarize: The 6 Steps of Your no!no! HAIR LITE Treatment

WARNING! Never use no!no! HAIR LITE on the genitals, nipples or face.

1. Remove the safety seal
2. Plug it in
3. Turn on the unit by pushing the on/off button on the top of the unit.
4. Glide in one direction in a slow, steady motion – always pulling no!no! HAIR LITE upward towards your body with the signal light facing you.
5. Work on one small section of the treatment area at a time, completing one before moving on to the next. Go over each area 3-4 times.
6. Complete the treatment by buffing the crystallized hair with the no!no! Buffer.

Thermicon™ Tips

To remove the Thermicon™ Tip, grasp it by both ends and gently pull it out of the socket. To replace it with a new tip, grasp it by both ends, line up the two small holes on the Thermicon™ Tip with the posts inside the Thermicon™ Tip socket and push it gently but firmly into place. Do not use force.

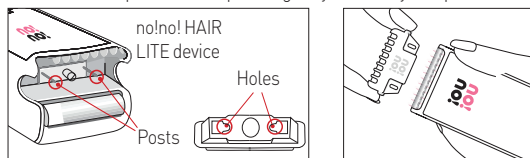


Figure 11: Installing and removing the Thermicon™ Tip

Check the Thermicon™ Tip to make sure it is inserted correctly. Make sure the Thermicon™ Tip is pushed all the way into the socket.



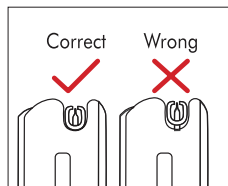


Figure 12

Always check the Thermicon™ Tip to make sure it is intact.

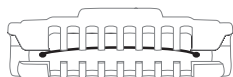

Bent thermodynamic wire

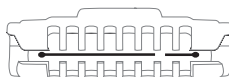

Torn thermodynamic wire

Figure 13

Maintenance

Do not use any liquids of any kind to clean the Thermicon™ Tip or rollers.

Thermicon™ Tips

It is important to keep the Thermicon™ Tip free of hair debris. Allowing hair to build up may cause it to burn out prematurely or malfunction. Check and clean the tip periodically during each treatment session.

To clean the Thermicon™ Tip:

1. Turn off your no!no! HAIR LITE
2. Remove the Thermicon™ Tip
3. Gently brush clean (Figure 14)

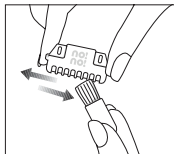


Figure 14

Rollers

Although not necessary to the performance of no!no! HAIR LITE, keeping the rollers clean is all part of proper maintenance and will help the tips last longer. Use the cleaning brush to remove excess hair and debris from the rollers.

Brush debris out and away from the socket in order to prevent any particles from falling into the socket.

If any hair or debris fall into the socket, this may cause no!no! HAIR LITE to not function properly. Gently blow any particles out of the socket before beginning the next treatment.

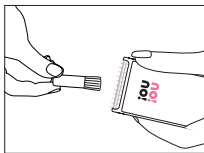


Figure 15

Troubleshooting

If no!no! HAIR LITE does not turn on:

- Make sure that the power cord is plugged firmly into the device and the electrical wall outlet.
- Make sure the On/Off button is switched to On.
- Make sure the Thermicon™ Tip is intact. Replace the Thermicon™ Tip if necessary.
- Replace the Thermicon™ Tip if the Thermodynamic Wire is severed
- Replace the Thermicon™ Tip if the Thermodynamic Wire is bent.

If the red light turns on:

- Replace the Thermicon™ Tip

If the device is not working properly:

- Contact your point of purchase for examination or repair.

The no!no! HAIR LITE warranty is valid for one year and only if the device has been used according to the instructions in this guide. Retain all receipts and original packaging.

Technical Information:

Technical Specifications:

Weight: 70gr

Physical Dimensions: 3.8 X 3.3 X 81mm

(1.5 X 1.3 X 3.2 in)

Mode of Operation: Continuous

Operating Conditions:

Temperature: 5° C-30° C

Humidity: 30%-80% RH

Rating: 100-240VAC 50-60 HZ 0.3A

Transport and Storage


Temperature: -20° C to 80° C

Humidity: 05-95% RH


Atmospheric pressure: 0.1-1.5atm

Electrical Safety:

Class II equipment – double insulation

( power supply only)

This appliance is not suitable for use in the presence of flammable anesthetic mixture with air or with Oxygen or Nitrous Oxide.

The symbol  on the label affixed to this appliance means "Attention, consult accompanying documents".

The Quality Management System complies with the Quality Management Standard ISO 13485.

no!no! HAIR LITE complies with IEC 60335 and with Low Voltage Directive (Electric Safety)
LVD: 2006/95/EC

The design and technical specifications are subject to change without notice.

Manufacturer:

Radiancy Israel Ltd.
5 Hanagar Street
45240 Hod Hasharon
Israel

EC Representative:

Obelis S.A Av.de Tervuren 34 Bte. 44
B-1040 Brussels, BELGIUM
Tel: 32 (0) 2 732 59 54
Fax: 32 (0) 2 732 60 03
GSM: 07545 46 60
Email: mail@obelis.net
Issued patents: 6825445, 7170034.
Pending applications: US2006/0011024



www.my-no-no.com