

SALON CELEBRITY™

LASER HAIR REMOVER

LAHR2-C

instructions

WARNINGS

- **Read all instructions thoroughly and retain for future reference.**

- This appliance can be used by children aged 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning the use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and maintenance shall not be made by children without supervision.
- Not suitable for Afro-Caribbean or black skin complexions. Not suitable for Grey Hair.
- Laser hair removal should not cause discomfort. If you feel any discomfort or soreness the setting may be too high or your skin tone (colour) is not suitable for use with this product.
- Power levels should be selected based on your skin type as defined in the 'Fitzpatrick Skin Phototypes' section.
- Do NOT treat the following areas: Tattoos, permanent make-up, pimples, warts, hair grown on black spots, skin disease, birth marks, burns, lips, nipples, genital or anal areas, eczema, eyelashes, eyebrows, near the eyes, nasal or ear hair.
- Do NOT treat the same hair/area repeatedly in any one session.
- Inbuilt safety features protect the eye being exposed to the laser, however, never aim the Wand at or near eyes.
- Do NOT use after drinking alcohol or if you are feeling unwell or fatigued.
- Do NOT use near flammable materials.
- Reflection Hazard: Please remove all jewellery in vicinity of the treatment area before use.
- When not in use, store the key responsibly in a safe place away from children.
- Misuse of the unit can cause harm.
- Do NOT override any safety features or open the unit.



the beauty specialists

CAUTIONS

- Do NOT dismantle or repair. If you have any problems using the product, contact our Customer Care Department.
- Misuse of the unit can cause harm.
- Do NOT override any safety features or open the unit or wand.
- Only use this product with the Mains Adapter supplied.

SPECIFICATIONS



- THIS PRODUCT MUST BE USED RESPONSIBLY.



- PLEASE READ THESE INSTRUCTIONS CAREFULLY BEFORE USE.

- REFER TO YOUR USER MANUAL FOR:

- WARNINGS

- SECURITY BUTTON SEQUENCE

- HOW DO LASERS WORK?

- SUITABLE SKIN TYPES

- PATCH TEST

- HAIR GROWTH CYCLE

- WHY IS A SERIES OF TREATMENTS REQUIRED?

- YOUR QUESTIONS ANSWERED

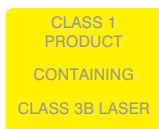
- TROUBLESHOOTING

808nm wavelength

Wearing laser safety goggles is unnecessary as the product has been designed with inbuilt safety features.



PATENTS PENDING



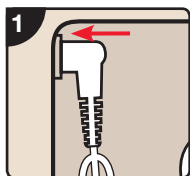
PATCH TEST

IMPORTANT: 24 hours before treatment it is very important to do a patch test.

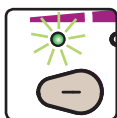
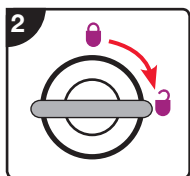
The patch test will confirm the maximum power level for your skin type. Each area to be treated must be checked as skin colour can vary.

Please refer to the 'PATCH TEST' section in the user manual for instructions.

- Best results are achieved on hairs that have been shaved a day or two before treatment or trimmed to 1mm–2mm in length.
- Remove all jewellery in the vicinity of the treatment area.
- Please ensure you have patch tested the area to be treated 24 hours prior to treatment.

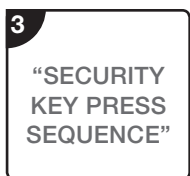
**STEP 1**

Plug Mains Adapter into socket on the Unit and into the Mains.

**STEP 2**

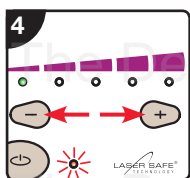
Use the Key to unlock the Keyswitch.

The left Power Level light will flash green.

**STEP 3**

Enter "SECURITY KEY PRESS SEQUENCE" within 60 seconds.

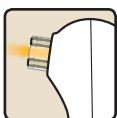
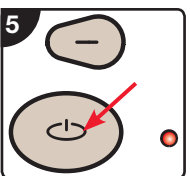
You will find the Security Key Press Sequence in the User Manual. Please ensure you read the User Manual before proceeding.

**STEP 4**

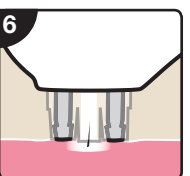
Use + and - buttons to set power level. Refer to User Manual for skin types.

If you feel any discomfort reduce the power level.

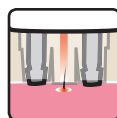
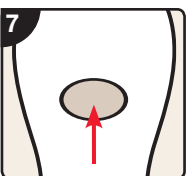
If the power level is adjusted the red Enable light will flash.

**STEP 5**

Press the Enable Button. The red Enable light will stay on and the amber light on the Wand will illuminate.

**STEP 6**

When placed on the skin the amber light will go OFF and the white 'ready' light will illuminate.

**STEP 7**

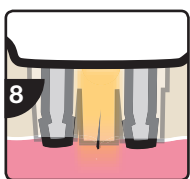
Press and HOLD the Laser Button to start the red laser treatment. You will hear a series of high pitch beeps to indicate the laser is operating. The laser appears as a small red light.

This must be aimed at the hair ROOT.

Note: On the lower power levels, the laser will pulse on and off for 4 seconds.

Note: You may feel a 'pin prick' sensation during treatment, however, most people do not feel anything.

Note: Sometimes you can smell faint burning or even see a very, small line of smoke as the hair is treated.

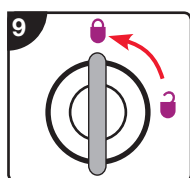


STEP 8

When the red laser light treatment ceases the amber light will come back on. You can now move onto the next hair.

Note: The Touch Pins must be removed briefly from the skin between treating each hair.

Note: Successfully treated hairs do not always fall out immediately, but are shed over time as they are released from the hair follicle, light rubbing can encourage this.



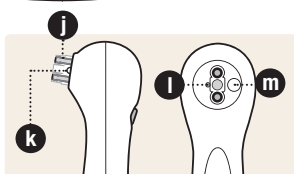
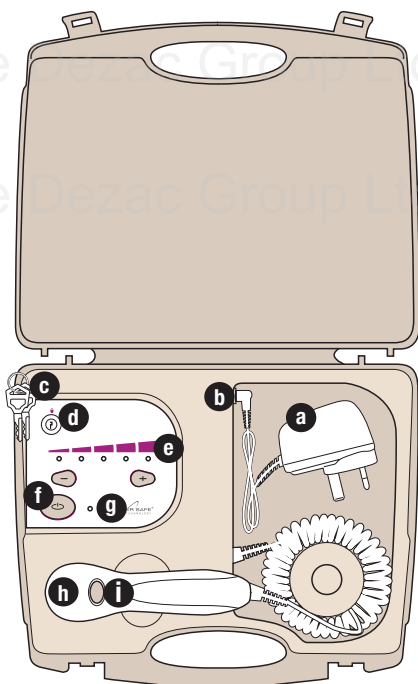
STEP 9

After use, turn the Keyswitch to lock the unit. Remove the Key to a safe place away from children.

Note: If no buttons are pressed after 20 minutes the Unit will enter energy saving mode. The orange light will flash and power level 1 will remain on. To restart turn key off and back on again and start at Step 3.

If you experience unexpected behaviour of this product then it is possible some external electromagnetic interference is responsible. You should remove this source and restart the product to restore full operation. At no time will the laser fire unintentionally.

COMPONENTS



- a. Mains Adapter
- b. Power Socket
- c. Keys
- d. Keyswitch
- e. Power Level (5 green lights)
- f. Enable Button
- g. Enable Light (red)

- h. Hand-held Wand
- i. Laser Button
- j. Laser Shroud and Touch Pins
- k. Laser Aperture
- l. White 'ready' Light
- m. Amber Light