


User Manual

Release 1.0

Electrochemical Dilatometer ECD-3-nano

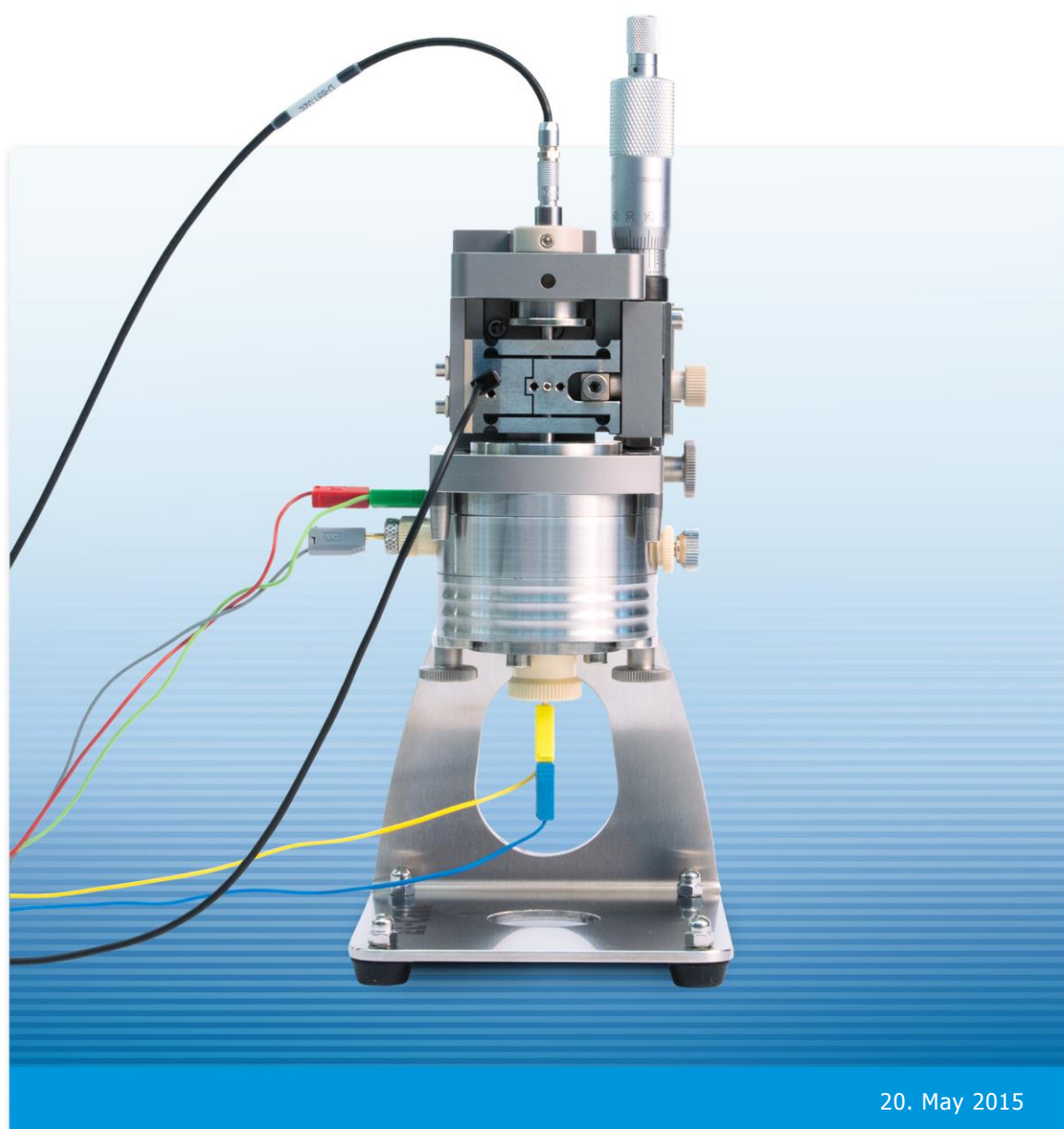

20. May 2015

The information in this manual has been carefully checked and believed to be accurate; however, no responsibility is assumed for inaccuracies.

EL-CELL GmbH maintains the right to make changes without further notice to products described in this manual to improve reliability, function, or design. EL-CELL GmbH does not assume any liability arising from the use or application of this product.

EL-CELL GmbH

Tempowerkring 8
21079 Hamburg - Germany

phone: +49 (0)40 790 12 733

fax: +49 (0)40 790 12 736

e-mail: info@el-cell.com

web: www.el-cell.com


Content

1	Product Description	4
2	Features	7
3	Safety Precautions	8
4	Unpacking	8
5	Start-Up	10
6	Assembling the cell inside the glove box	14
7	Further assembling outside the glove box	18
8	EC-Link Software Installation	22
9	Calibration and Settings	22
10	Recording the Displacement Signal with an External Potentiostat	23
11	Using the Reference Electrode	24
12	Dilatometer Disassembly and Cleaning	25
13	Care Instructions	26
14	Consumables	27
15	Technical Support	28
16	Warranty	28
 Appendix:		
	Components Sensor Unit	29
	Components Cell Body	30
	Connector and Cable Pin-out	31

1 Product Description

The ECD-3-nano electrochemical dilatometer is dedicated to the measurement of charge-induced strain (expansion and shrinkage) of electrodes down to the nanometer range. The ECD-3-nano has been particularly developed for the investigation of Li-ion battery and other insertion-type electrodes. It may, however, also be used for many other electrochemical systems utilizing aprotic organic electrolyte solutions. The electrode materials used can either be bound film or single crystals/grains (e.g. graphite flakes). The maximum sample size is 10 mm x 1 mm (diameter x thickness).

The heart of the ECD-3-nano is an electrochemical cell, hermetically sealed against ambient atmosphere. The two electrodes inside are separated by a stiff glass frit which is fixed in position. The upper working electrode (**WE**) is sealed by means of a thin metal foil, through which any charge-induced thickness change is transmitted towards the sensor/load unit above. This working principle allows determining the height change

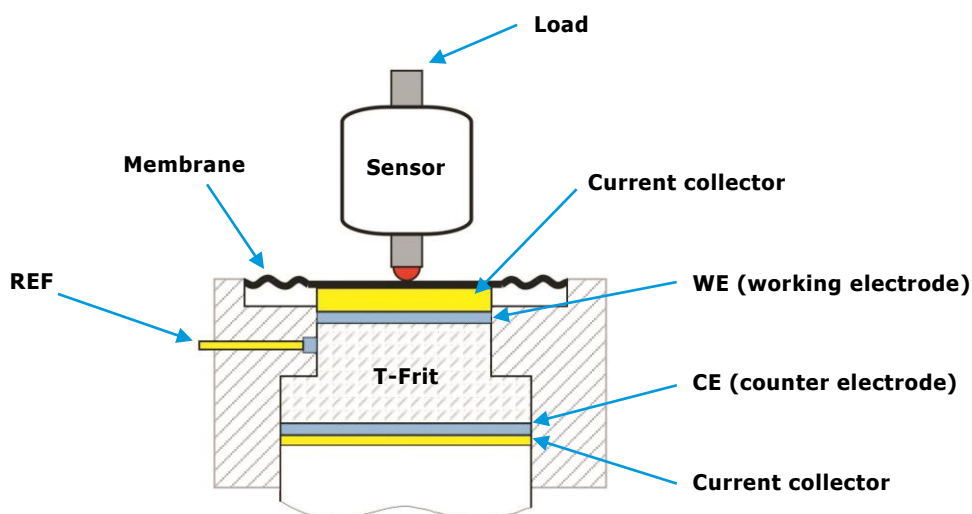


of the working electrode without any interference from that of the counter electrode (**CE**).

A high-resolution capacitive displacement transducer detects dimensional changes of the WE ranging from a few nanometers up to 250 micrometers during one and the same experiment that may last between a few minutes to many days.


The ECD-3-nano features an integrated USB data logger for recording the electrode displacement, temperature, cell potentials and current. Analog outputs of displacement and temperature are provided for integration with external instruments.

For best accuracy and drift stability, the dilatometer is to be operated inside a temperature controlled chamber.

Basic structure of the ECD-3-nano:

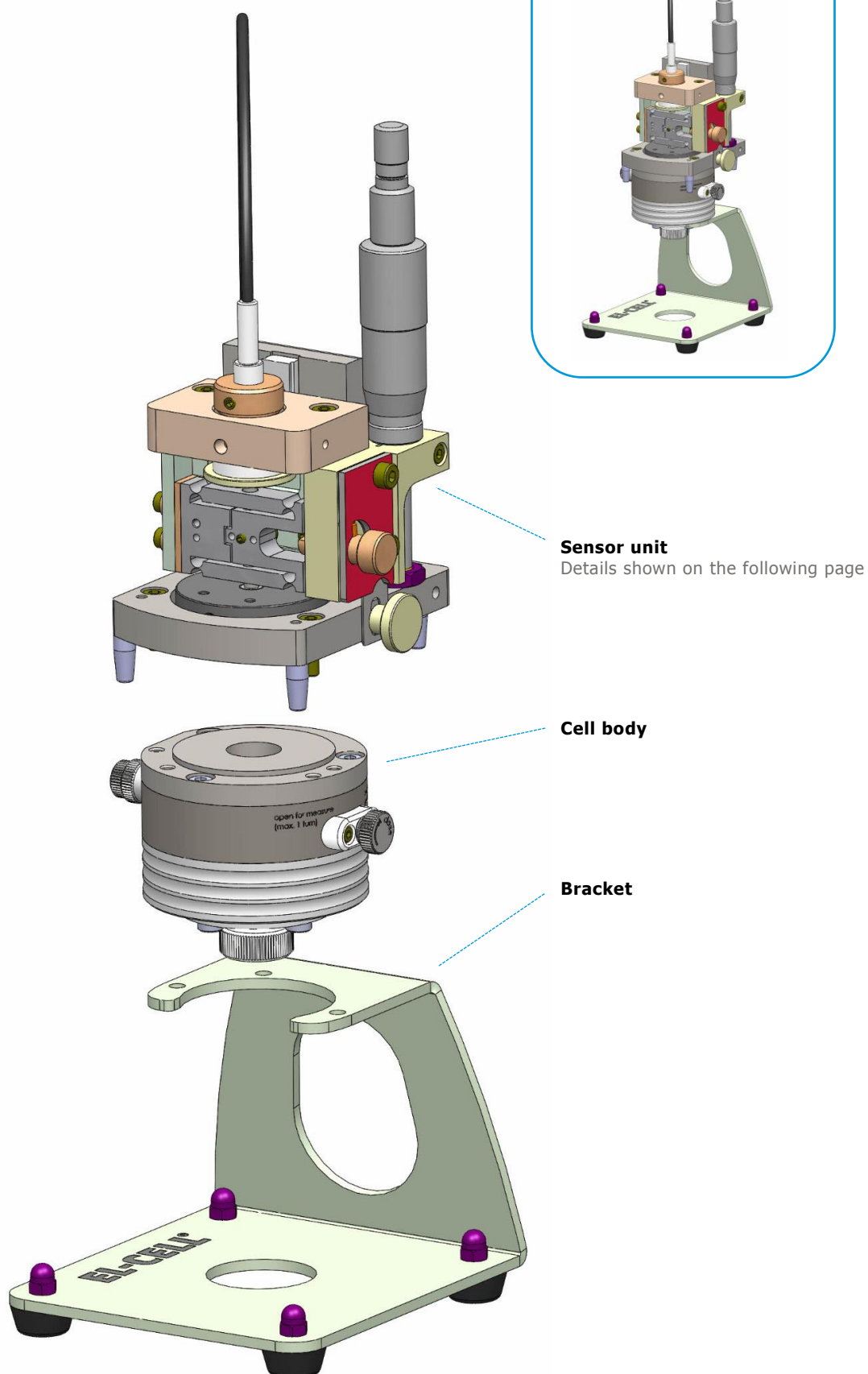

Cut drawing of the ECD-3-nano:

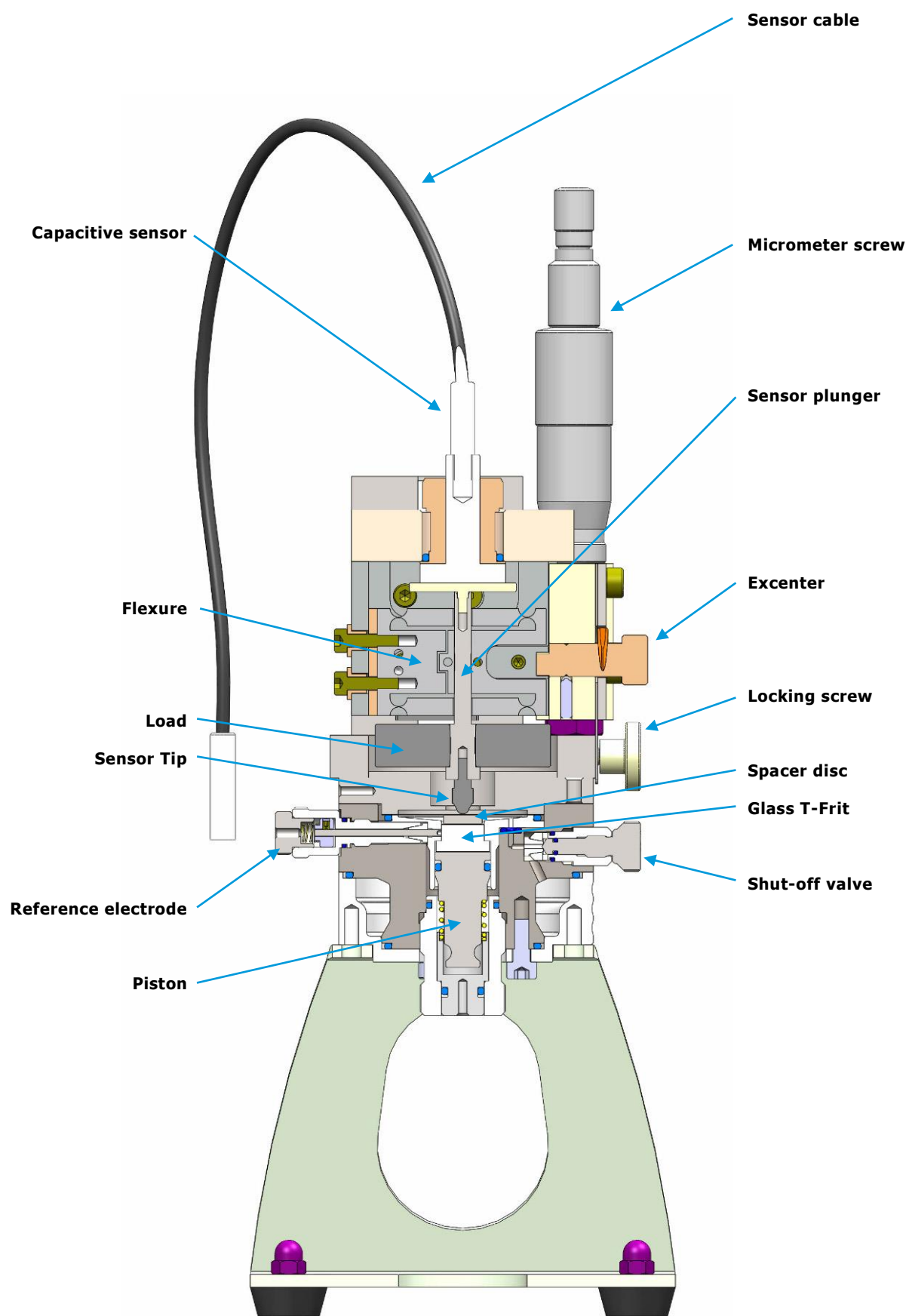


2 Features

The ECD-3-nano is an electrochemical dilatometer for measuring changes of thickness of the working electrode of a battery test cell. The main features of the ECD-3-nano are briefly described in the following:

- High resolution capacitive sensor system with <5 nm resolution, drift stability of <20 nm/hour (sample-free instrument at constant temperature), and 250 μm full range.
- Conditioning electronics with analog output signals (-10 to 10 V) for displacement and temperature.
- Integrated USB data logger for the recording of displacement, temperature, cell potentials and current.
- 3-electrode electrochemical cell
- Sample (working electrode): bound electrode film or single crystal/grain; max. sample size 10 mm x 1 mm (diameter x thickness)
- Load on working electrode: approx. 1 N
- Electrolyte volume: approx. 1 ml
- Materials in contact with electrolyte: PEEK, borosilicate glass, EPDM rubber, stainless steel 316L for aprotic, gold for aqueous electrolytes
- **Operating temperature range:** Cell and sensor: -20 to +70 °C; Conditioning electronics and data logger: 0 to + 40 °C

Dimensions of the ECD-3-nano (mm):

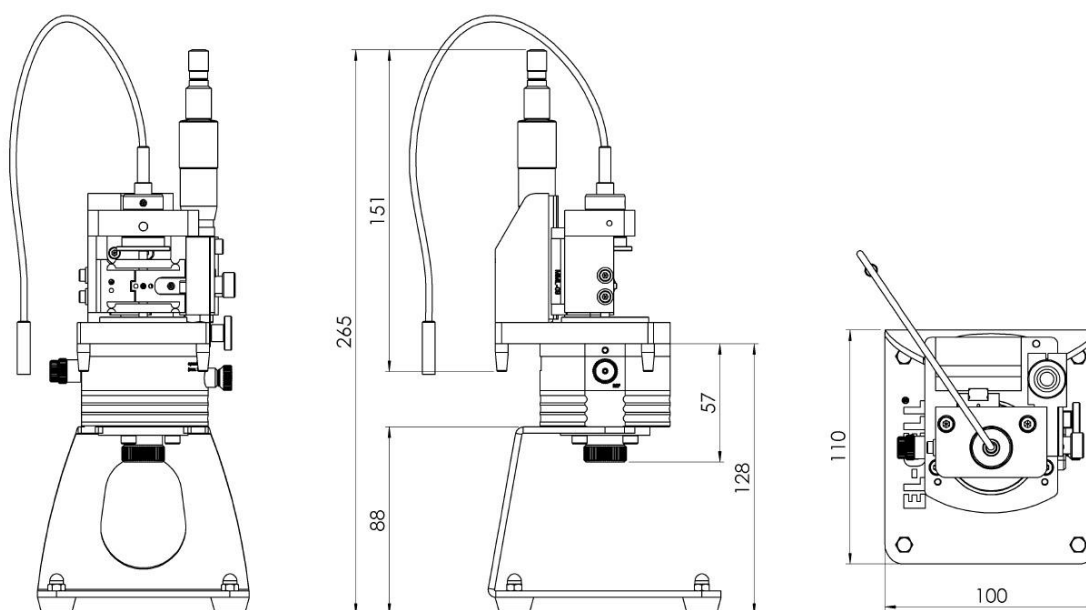

3 Safety Precautions


Use proper safety precautions when using hazardous electrolytes. Wear protective glasses and gloves to protect you against electrolyte that may accidentally spill out of the instrument during filling, operation and disassembly.

4 Unpacking

Check the contents of the packages against the list given below to verify that you have received all of the required components. Contact EL-CELL, if anything is missing or damaged. **NOTE:** Damaged shipments must remain within the original packaging for freight company inspection.

List of Components:

1. ECD-3-nano dilatometer [ECD3-00-0026-A](#), assembled
2. Box ECD-3-nano [ECE1-00-0006-F](#), assembled
3. Sensor-GND cable [ECE1-00-0041-A](#)
4. Sensor cable (PISeca) [SEN9023](#)
5. ECD cell cable [ECE1-00-0033-E](#)
6. Power supply SPU 45E-303 [ELT9207](#)
7. Cable connector NC5FRX (5 pole, female) [ELT9207](#)
8. USB cable typ A/B (2.0 m) [ELT9167](#)

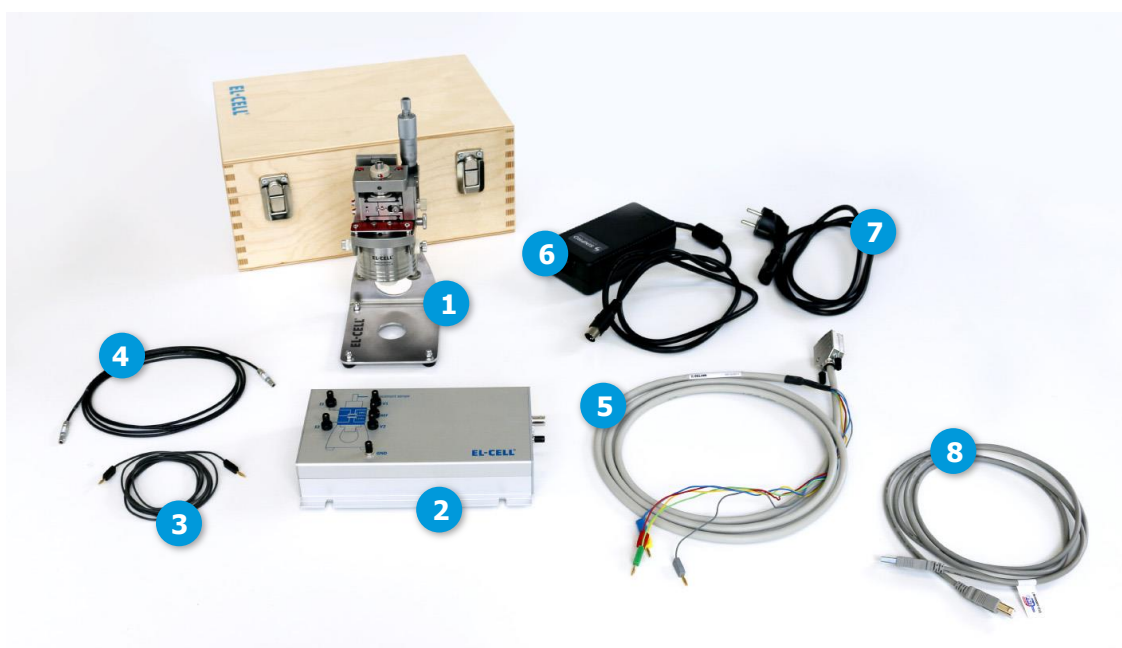

Accessories Kit: ECD3-00-0400-A

1. 4 x O-Ring 33.05 mm x 1.78 mm [DIC9005](#)
2. 2 x O-Ring 9.75 mm x 1.78 mm [DIC9006](#)
3. 2 x O-Ring 2 mm x 1.5 mm [DIC9025](#)
4. 2 x O-Ring 50.5 mm x 1.78 mm [DIC9026](#)
5. Membrane (aprotic) 1.4404 [ECC1-00-0019-D](#)
6. Spacer disc (set) 2.1 - 2.3 [ECC1-01-0012-F](#)
7. Demonstration kit (5 x activated carbon electrode foil w/5% PTFE Binder, 10 mm) [ECD1-00-0900-A](#)
8. Filling tube [ECD3-01-0001-A](#)
9. CD containing EC-Link data logger software [ECE1-00-0052-A](#)
10. Tweezers [WZG9001](#)
11. Spherical allen screw driver 3 mm [WZG9002](#)
12. Allen screw driver 2.5 mm [WZG9003](#)
13. Vacuum tweezers [WZG9004](#)
14. Set allen wrench [ECC1-01-0028-A](#)




5 Start-Up

If you are using the ECD-3-nano the first time, please follow the below procedure starting at step 1.



After unpacking the ECD-3-nano, remove the transport lock from the sensor unit.

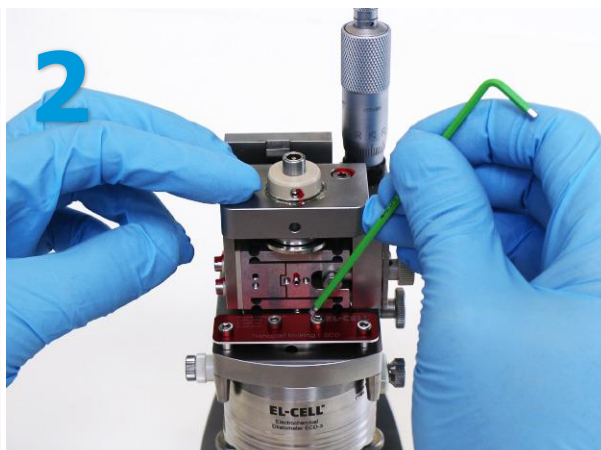

Remove the two inner screws first, then the two outer screws. After that, screw off the sensor unit from the bracket.



Screw off the cell body from the bracket.




Now unscrew the spring load from the bottom of the cell body.




Unscrew the reference electrode.




Remove the cover flange by loosening the three screws.

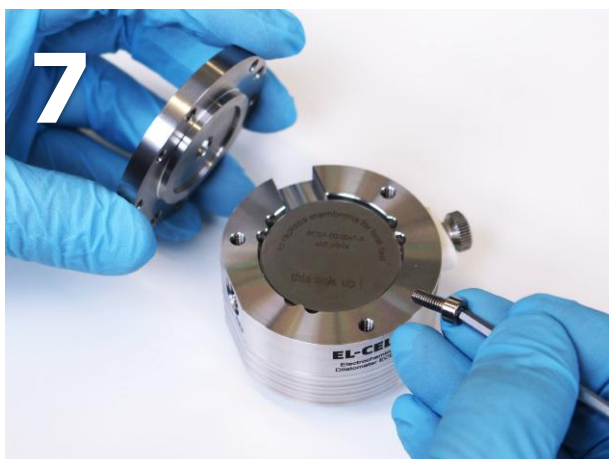

There is a stiff plate under the cover flange, which must later on be replaced by the provided flexible membrane.

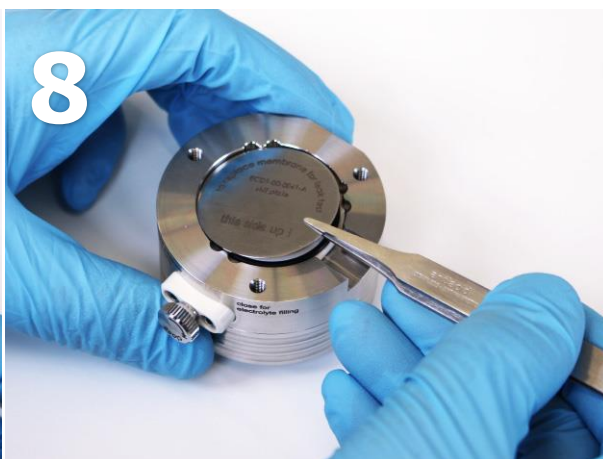

Remove the stiff plate from the cell body. Later the stiff plate can be used for checking the cell for leaks if necessary.

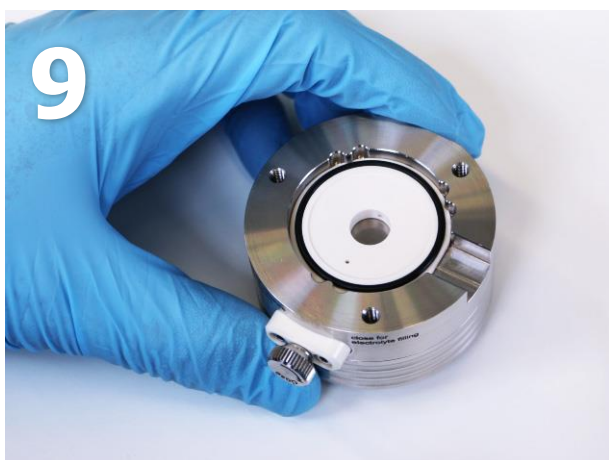

Now the frit flange with the O-Ring and the piston in the middle are visible.



Push the frit flange out of the cell body.



11

Make sure not to loose the little O-Ring seal.



12

Pull the piston out of the frit flange. Remove the T-Frit afterwards.



13

Remove the dead volume cover from the ECD-3 base body by unscrewing the three screws at the cell bottom.



14

Remove both O-Rings.


All the below shown parts need to be dried before they can be moved into the glove box for assembly. Recommended drying conditions: 80°C, <0.01 mbar, 12 hours.



1. Membrane (aprotic) 1.4404
2. Spacer disc (height depends on the height of your electrode)
3. T-Frit
4. Cover flange with three screws
5. Dead volume cover with three screws
6. ECD-3 base body
7. Frit flange
8. Spring load
9. Reference electrode
10. Piston
11. O-Ring 50.5 x 1.78 mm
12. 2 x O-Ring 33.05 x 1.78 mm
13. O-Ring 2 x 1.5 mm

6 Assembling the cell inside the glove box

After moving the different parts of the disassembled cell body into the glove box, follow the steps below. Protect yourself and handle the chemicals with care.

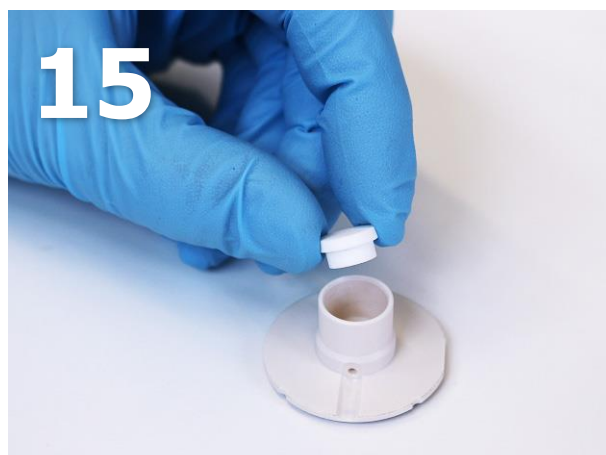

Inside the glove box: Insert the T-frit with the smaller side pointing downwards into the frit flange.



Inside the glove box: Put the lithium counter electrode (10 mm diameter) on top of the T-frit into the frit flange.




Inside the glove box: Attach the counter piston from below.



Inside the glove box: Insert the two big O-rings and assemble the dead volume to the cell base again.



Inside the glove box: Put this assembly into the cell base body. Make sure to have the two grooves properly aligned. Don't forget to insert the little O-Ring seal!

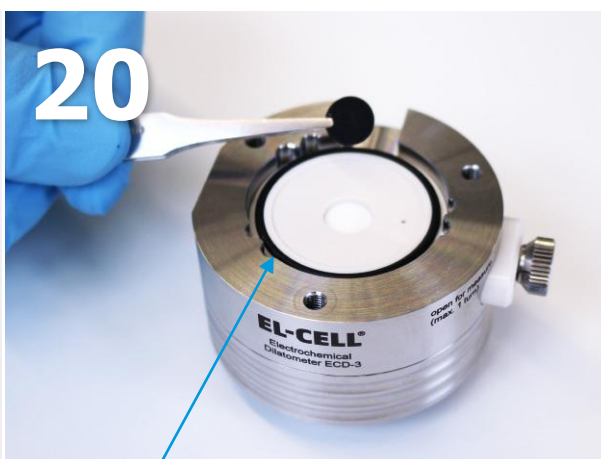

Inside the glove box: Insert the O-ring seal (DIC9005, see arrow). Then place the working electrode with the active side down on top of the T-Frit.




Inside the glove box: Put one of the provided spacer discs on top of the electrode. Use the 2.3 mm thick spacer disc for samples of $<100\ \mu\text{m}$ thickness.



Inside the glove box: Then put the membrane on top.



Inside the glove box: Attach the cover flange.

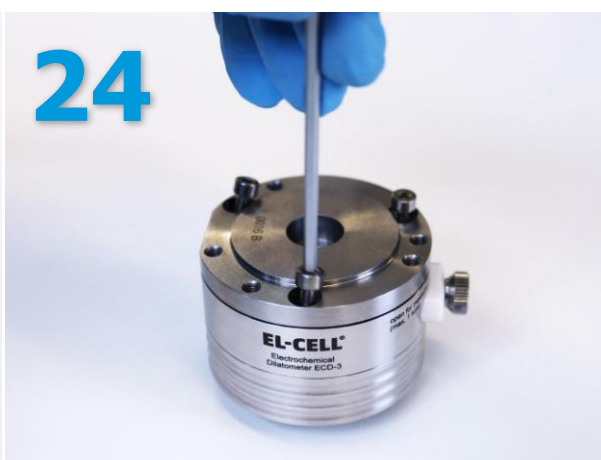

Inside the glove box: Close the cell body by tightening the three screws.



Inside the glove box: Now screw in the spring load into the cell base.



Inside the glove box: Before filling the cell body, close the shut-off valve clockwise.



Inside the glove box: Load the syringe with approx. 1.5 ml of electrolyte and connect the syringe to the cell body.



Inside the glove box: Pull back the syringe piston in order to evacuate the cell. Hold the vacuum for a few seconds. Then release the piston.



Inside the glove box: The electrolyte will be sucked into the cell by the vacuum applied. Never push the syringe piston! Then remove the fill line and syringe.



Inside the glove box: Pick up some lithium with the reference pin. Make sure that the hole of the reference pin is completely filled with lithium metal.



31

Inside the glove box: Lithium must not come into contact with the PTFE ferrule (see arrow)!



32

Inside the glove box: Attach the reference pin to the cell body.



33

Inside the glove box: Make also sure that there is no lithium outside the pin before attaching it to the cell body.




34

Inside the glove box: The cell is now assembled and hermetically sealed, so that you can take it outside the glove box.

7 Further assembling outside the glove box




Hook the assembly into the bracket and fasten it with the two knurled screws.

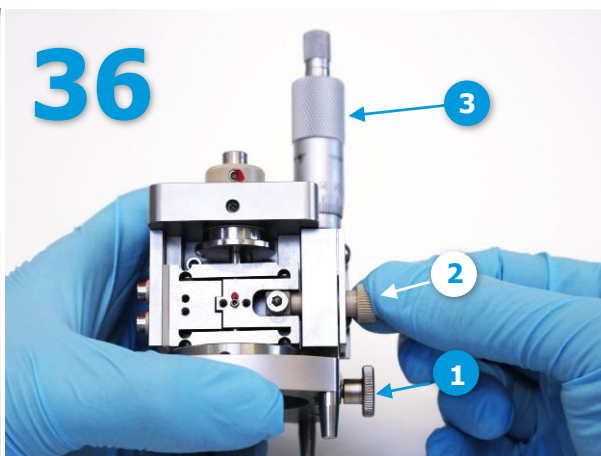

Release the locking screw (1) and the excenter (2) of the sensor unit. Then move the sensor tip into the upmost position by turning the micrometer screw (3) clockwise.




Attach the sensor unit onto the dilatometer cell.

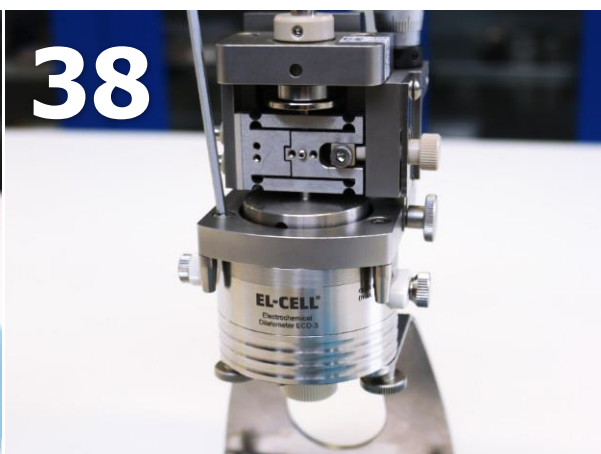

Fasten the screws to fix the sensor unit on top of the cell.


Now connect all cables as shown in the pictures

1. ECD cell cable
2. Sensor GND cable
3. Sensor cable
4. Power cord (90 to 240V AC input)
5. DC cable (+15V, -15V, +5V DC output)
6. USB 2.0 cable
7. Optional analog output cable

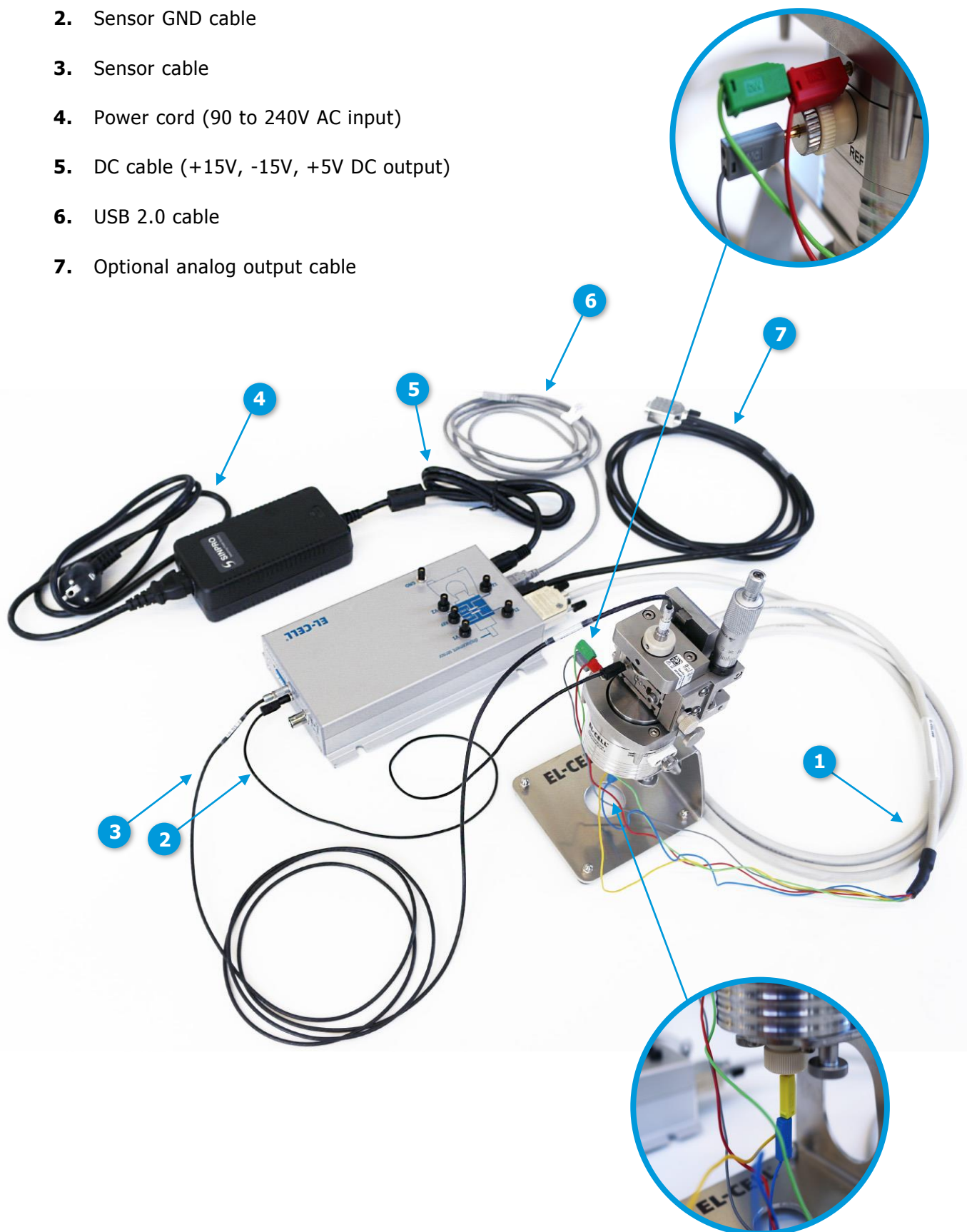




Release the locking screw (1) and the excenter (2) if not already done.

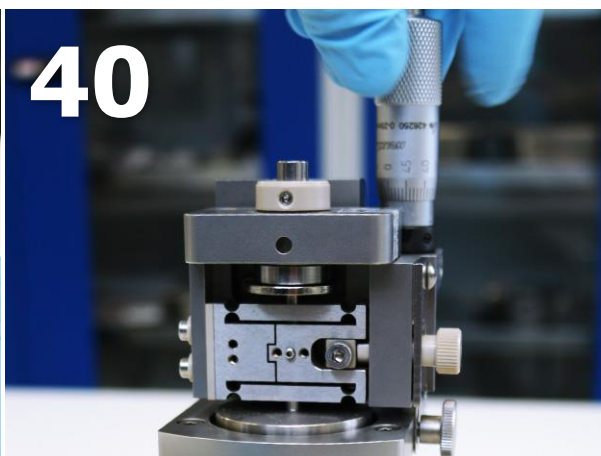


Adjust the sensor position by turning the micrometer screw counter clockwise.

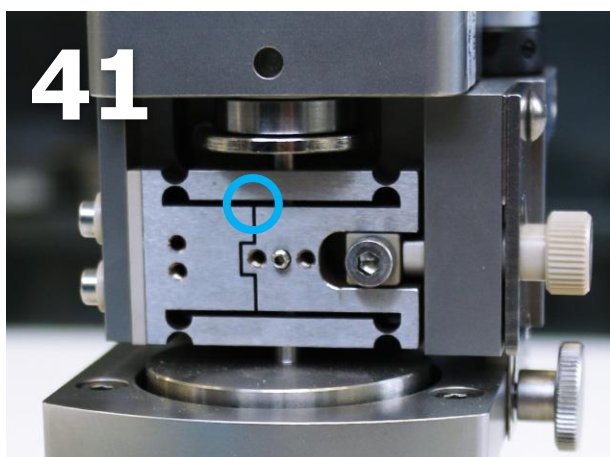

The sensor position can be viewed at the flexure guide ..



.. and at the bar graph indicator of the controller box. Adjust the sensor position around mid-gap.

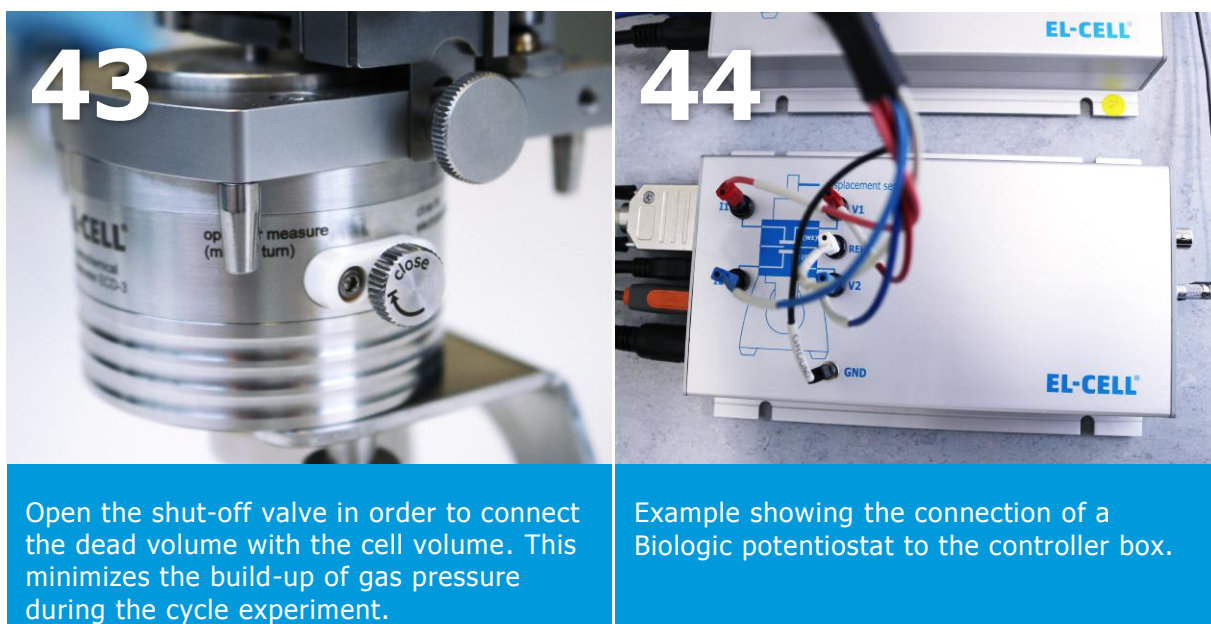

Open the shut-off valve in order to connect the dead volume with the cell volume. This minimizes the build-up of gas pressure during the cycle experiment.

Example showing the connection of a Biologic potentiostat to the controller box.

Finally, connect your potentiostat or battery tester to the 4 mm jacks on the front panel of the controller box. The rightmost column in the table below refers to the terminology used for the lead connections of Biologic potentiostats (MPG-2, SP, VSP and VMP series). <http://www.bio-logic.info/electrochemistry-ec-lab/instruments/>

Controller Box	Potentiostat	Biologic Potentiostat, VSP, VMP3, etc.
I1	WE Current	WE
V1	WE Sense	Ref1
REF	Reference	Ref2
V2	CE Sense (if available)	Ref3
I2	CE	CE
GND	GND (if available)	GND

Before starting the electrochemical cycle we recommend holding the cell at constant potential (or open-circuit) for several hours to allow for baseline stabilization. The initial rest period helps to discern charging induced dimensional changes from the initial creeping.

NOTE: All materials display a more or less pronounced creeping. They tend to shrink when applying a load, and to swell when removing this load. A mayor contribution to the initial creeping seen right after cell assembly is to be assigned to the construction materials of the dilatometer. Creeping of the working electrode is induced each time the mechanical properties of the working electrode are altered by charging. Therefore, each charge induced height change is followed by some creeping. The charge induced creeping effects are real and not artefacts of the measurement.

8 EC-Link Software Installation

In order to record the displacement signal together with the cell voltage, cell current, electrode potential and temperature, the software of the integrated data logger needs to be installed on a Windows® PC.

- You must be logged into an account with Administrator privileges.
- Save your work and close down all active programs.
- On the installation CD, run [X:\Driver_CDM20814_Setup](#) (where X refers to the CD drive). This will install the FTDI driver required to establish the USB connection with the data logger.
- On the installation CD, run [X:\setup](#). This will install the data logger software. Follow any instructions that may appear on your screen.
- Once installation is finished plug in the provided USB cable into both the host PC and the ECD-3-nano controller box.
- Launch the data logger software if not already done.
- After a few seconds, the data logger software should report a valid connection and you are ready to start the measurement.

Additional information on the EC-Link software can be found in a separate manual (<http://el-cell.com/downloads/downloads-manuals>).

9 Calibration and Settings

Calibration of the instrument has been carried out at the factory. The corresponding settings of the EC-Link software are stored in the file [Settings ECD-nano-3 \[Device-ID\].txt](#) in the installation directory on the local hard drive and on the installation CD. If the default settings have been changed for any reason, the original settings can be restored by copying [Settings ECD-nano-3\[Device-ID\].txt](#) from the installation CD into the directory C:\ProgramData\EC-LINK. The settings affect only the data logger readings.

The DIP switches at the controller box affect both the displacement readings of the USB data logger and the analog output signal. Leave the switches in their default position as shown in the screenshot below.



10 Recording the Displacement Signal with an External Potentiostat

Many of today's battery testers and potentiostats provide additional analog inputs that may be used to record sensor signals along with cell current and potential.

In the following, the combination of the ECD-3-nano with a Biologic potentiostat (MPG-2, SP, VSP and VMP series) is described as an example. The Biologic potentiostats feature two analog inputs that are used here to record both displacement and temperature.

1. Connect the 9-pin Sub-D connector of the optional analog output cable to the analog input of the respective VMP3 channel.
2. In the Biologic EC-Lab software, load the experiment settings [ECD-3-nano.mps](#) provided on the ECD-3-nano documentation CD. The settings are shown in the External Devices dialog (see screenshot below; actual settings may differ). Adapt the [Parameter Settings](#) of the charge/ discharge protocol to your particular experiment, if necessary.

The screenshot displays the 'External Devices' dialog box in the Biologic EC-Lab software. It contains two sections for configuring analog inputs:

- Analog IN 1:** The 'Convert E/V' checkbox is checked. The 'to' dropdown is set to 'Displacement/ μm '. Below, the 'with' section shows two rows: the first row has '10' in the input field, 'V =', '-125' in the input field, and ' μm (max)'; the second row has '-10' in the input field, 'V =', '125' in the input field, and ' μm (min)'.
- Analog IN 2:** The 'Convert E/V' checkbox is checked. The 'to' dropdown is set to 'T/ $^{\circ}\text{C}$ '. Below, the 'with' section shows two rows: the first row has '0' in the input field, 'V =', '0' in the input field, and ' $^{\circ}\text{C}$ (max)'; the second row has '10' in the input field, 'V =', '80' in the input field, and ' $^{\circ}\text{C}$ (min)'.


11 Using the Reference Electrode

The reference electrode assembly is comprised of the reference electrode, the set collar attached to the pin by means of a set screw, the fitting, the spring, and the hollow screw, cf. the sketch below.

The hollow screw serves to apply the spring pressure on the set collar, thereby gently pushing the reference pin against the glass frit. The blind bore on the tip of the reference pin is intended for taking up the reference electrode material. For most lithium ion chemistries the reference material may be a small piece of lithium metal picked up by the reference pin. For other aprotic electrolytes, and also for some aqueous systems, a piece of PTFE bound activated carbon may serve as the (pseudo) reference material. The optional gold reference pin is recommended for use in aqueous electrolytes.

NOTE: Do not use the gold reference pin in combination with lithium metal as the reference material.

Components of the reference electrode:

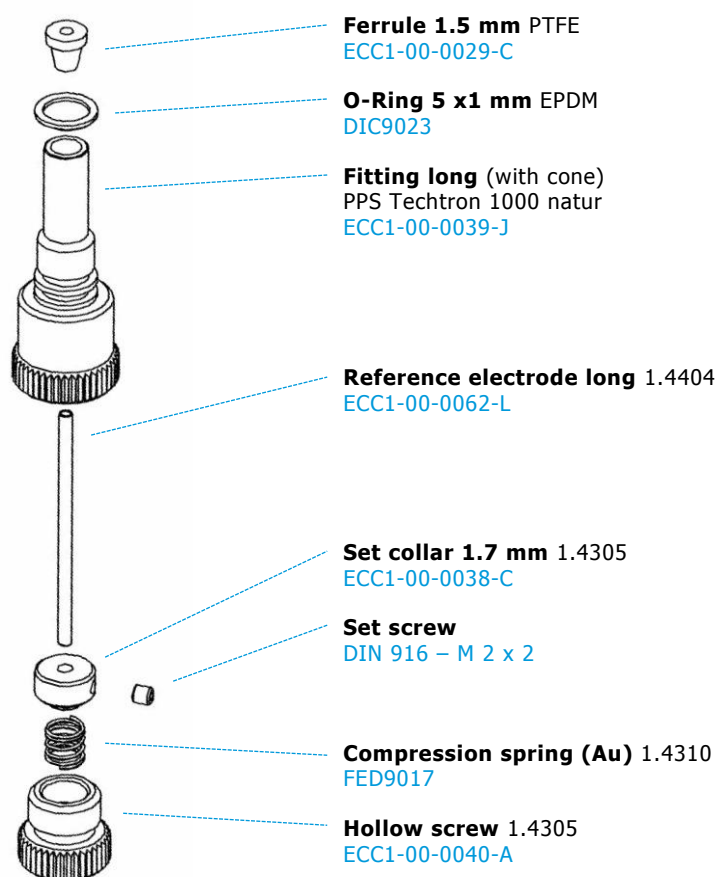

12 Dilatometer Disassembly and Cleaning

When disassembling the dilatometer cell, wear protective gloves and glasses.

Collect parts that have been in contact with electrolyte on a separate tray for subsequent cleaning.

1. Disconnect all cables from the dilatometer cell and the sensor unit.
2. Remove the dilatometer cell from the temperature chamber.
3. Detach the sensor unit from the dilatometer cell.
4. Detach the cell from the bracket.

Then disassemble the cell body as described in chapter 5, starting at step 4. It is very important to unscrew the spring load first. After unscrewing the spring load you can detach the reference electrode.

In contrast to the initial start-up of the dilatometer, also the shut-off valve must be disassembled after testing.



Important: Dissassemble and clean the shut-off valve.

Clean all wetted parts right after disassembly. Ultrasonic cleaning with water and/or detergent wash is recommended. Valves and tubing may clog if not properly purged with water or other solvent. After cleaning, dry all parts in vacuum at 80°C overnight.

For further tests follow the instructions beginning at step 15, chapter 6.

13 Care Instructions

Maintain the O-Rings regularly with silicone grease. We recommend to use new O-Rings for each measurement in order to achieve best results.



Take some tiny amount of silicone grease on a sponge or cloth.



Wipe the seal at the valve stem with the sponge..



.. and wipe the O-ring seals.

14 Consumables

Cell Body:

- T-Frit 10/12.5 [ECC1-00-0041-B](#)
- Membrane (aprotic) 1.4404 [ECC1-00-0019-D](#)
- O-Ring 33.05 x 1.78 mm [DIC9005](#)
- O-Ring 2 x 1.5 mm [DIC9025](#)
- O-Ring 50.5 x 1.78 mm [DIC9026](#)
- O-Ring 9.75 x 1.78 mm [DIC9006](#)
- Accessories kit ECD-3 [ECD-3-00-0400-A](#)

Sensor Unit:

- Socket screw [DIN-912 M4 x 12](#)

15 Technical Support

Technical support for this product is exclusively provided by EL-CELL GmbH.

EL-CELL GmbH

Tempowerkring 8

21079 Hamburg - Germany

phone: +49 (0)40 790 12 733

fax: +49 (0)40 790 12 736

e-mail: info@el-cell.com

web: www.el-cell.com

16 Warranty

For a period of one year from the date of shipment, EL-CELL GmbH (hereinafter Seller) warrants the goods to be free from defect in material and workmanship to the original purchaser. During the warranty period, Seller agrees to repair or replace defective and/or nonconforming goods or parts without charge for material or labor, or, at the Seller's option, demand return of the goods and tender repayment of the price. Buyer's exclusive remedy is repair or replacement of defective and nonconforming goods, or, at Seller's option, the repayment of the price.

Seller excludes and disclaims any liability for lost profits, personal injury, interruption of service, or for consequential incidental or special damages arising out of, resulting from, or relating in any manner to these goods.



This Limited Warranty does not cover defects, damage, or nonconformity resulting from abuse, misuse, neglect, lack of reasonable care, modification, or the attachment of improper devices to the goods. This Limited Warranty does not cover expendable items. This warranty is void when repairs are performed by a non-authorized person or service center. At Seller's option, repairs or replacements will be made on site or at the factory. If repairs or replacements are to be made at the factory, Buyer shall return the goods prepaid and bear all the risks of loss until delivered to the factory. If Seller returns the goods, they will be delivered prepaid and Seller will bear all risks of loss until delivery to Buyer. Buyer and Seller agree that this Limited Warranty shall be governed by and construed in accordance with the laws of Germany.

The warranties contained in this agreement are in lieu of all other warranties expressed or implied, including the warranties of merchantability and fitness for a particular purpose.

This Limited Warranty supersedes all prior proposals or representations oral or written and constitutes the entire understanding regarding the warranties made by Seller to Buyer. This Limited Warranty may not be expanded or modified except in writing signed by the parties hereto.



Components Sensor Unit

There are no further spare parts available for the sensor unit. For repair, please contact EL-CELL.

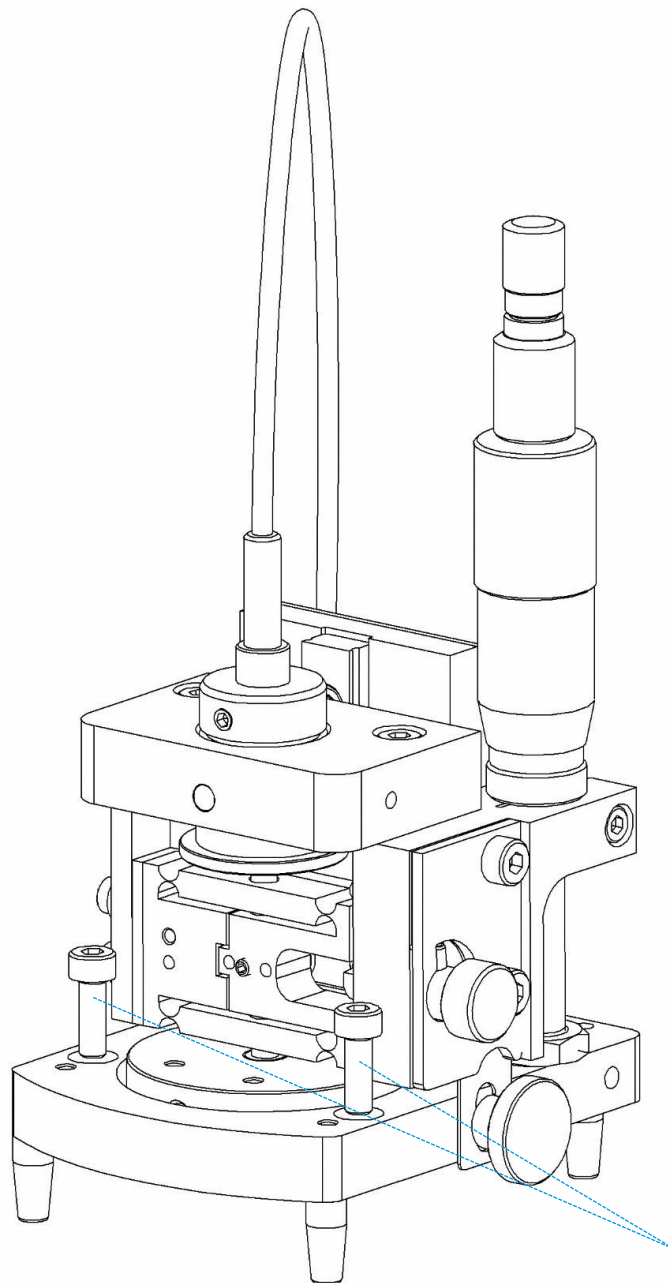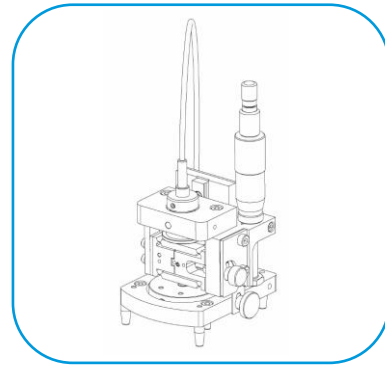

Socket screw
DIN-912 M4x12

Components Cell Body

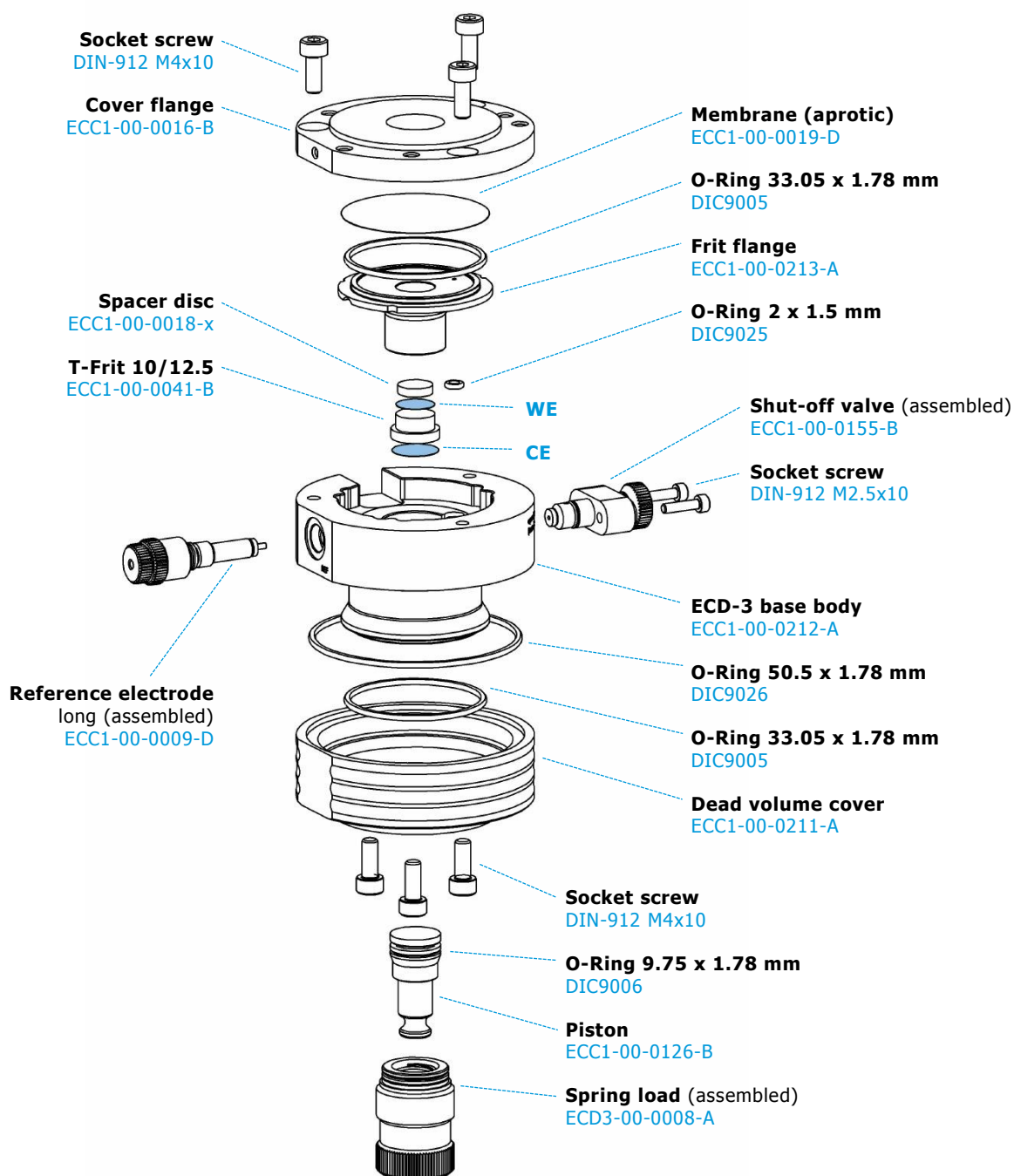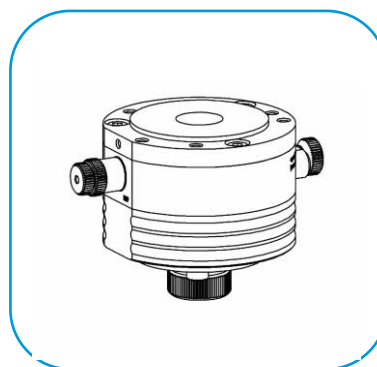

Connector and Cable Pin-out

Cell Cable (4 x 2 x 0.25 mm², TP, shielded):

ECE1-00-0033-E

One end of the cable is terminated with a Sub-D HD M15 connector (to box); the other end is terminated with 2 mm banana connectors. A Pt100 sensor is located beneath the black shrink tube at the end of the cable pointing to the dilatometer. The cable shield is connected to GND.

Pin #	Signal	Cable Color	Color of 2 mm connector
1	V1	Red	Red
2	V2	Blue	Blue
3	-	-	-
4	REF	Grey	Grey
5	I2	Yellow	Yellow
6	-	-	-
7	-	-	-
8	-	-	-
9	-	-	-
10	I1	Green	Green
11	Pt100(1)	Brown	-
12	Pt100(2)	White	-
13	-	-	-
14	-	-	-
15	-	-	-

Biologic Auxiliary Cable (2 x 2 x 0.14 mm², TP, shielded):

ECE1-00-0039-B

Both connector housings are tied to the cable shield. The cable shield is connected to GND.

IEEE 1394 to Box			Sub-D M9 to Biologic AUX Input	
Pin #	Signal	Cable Color	Pin #	Signal
1				
2	GND	Black	7	GND
3				
4	Temperature	Blue	6	Analog IN2
5				
6	Displacement	Green	1	Analog IN1