

**5in1 Mesoporation+LED Photon+Electroportion+EMS+RF**

**Radio Frequency Skin Lifting**

**SR-265**

**Use Manual**



## CATALOG

<b>Preface</b>	-----	3
<b>I Product Introduction</b>	-----	4
<b>II Features</b>	-----	5
<b>III Functions</b>	-----	6
<b>IV Specification</b>	-----	7
<b>V Package Listing</b>	-----	7
<b>VI How To Use</b>	-----	7
<b>VII Operation Steps</b>	-----	8
<b>VIII Forbidden Matter</b>	-----	12
<b>IX Daily Inspection And Maintenance</b>	-----	12

## **Preface**

Dear User:

Thank you for purchasing !Be sure to read this instruction manual carefully before use so that you can make full and correct use of your massager.Please keep this manual available for future reference after you have read it through.

**This product does not belong to medical,please don't use for disease treatment**

## **I Product Introduction**

This SR-265 skin care device is good for body & face daily therapy. It works on skin rejuvenation, skin toning, wrinkle removal, anti-age, skin renewal, anti-allergy, pore size reduction, increase skin elasticity, and improve skin absorption of cosmetic products, and so on. It uses technologies of electroporation, mesoporation, UV LED, Electronic Muscle Stimulation and RF.

SR-265 contains 2 pairs of bipolar electrical tips and 16 pieces of LED lamps. People feel a little bit warm, comfortable and relax during the treatment. No pain, no harm. Skin product (e.g.essence) should be applied. And the SR-265 would help skin to absorb the products more efficiently. Daily treatment is available. About 10 minutes treatment everyday can keep sharp body, smooth skin and young face.



### **a.EP( Electroporation)**

Introducing beauty essence into your skin

This beauty machine is used for collagen injection without a needle. Macromolecular components such as collagen and hyaluronic acid cannot easily penetrate the deep layer of the skin from the surface. But with our special technology called

"Electroporation" you can achieve a high level of effectiveness without any pain.

### **b.MP( Mesoporation )**

It helps essence into your skin deeply

Mesoporation is a combination of mesotherapy and electroporation. In Europe it is widely used as mesotherapy without a needle. It achieves a high level of penetration into the skin. It is reported that Mesoporation's penetration degree is about 10 times that of conventional therapy in a 1mm depth of skin and 27 times in 4mm. An amazing whitening and anti-aging effect is possible.

### **c.LED( Light Emitting Diode )**

For treatment of acne lesions

It works efficiently without damaging the skin by using an especially sensitive LED light developed exclusively for beauty equipment. It uses a blue LED with a wavelength of 415nm (+-10nm).

### **d.EMS( Electrical Muscle Stimulation )**

For lifting skin, weight loss, cellulite treatment

Four elements are used to ensure that the electrical current works randomly and stimulates the muscle tissue very gently causing a secondary muscular movement that preserves the skin's elasticity and suppleness.

### **e. RF( Radio Frequency )**

For detoxification, lymphatic drainage, improve blood circulation

Electrical energy tends to follow the shortest distance between two electrodes. When electricity runs through living tissue the high energy flow generates heat as a characteristic of the electrical energy and the resistance to electricity by the living tissue. This is another feature of Dr. Arrivo that adapts electricity to be gentle on the skin

## **II Features**


- a. Combined with 5 technologies with Electroporation, Mesoporation, LED, Electrical Muscle Stimulation & RF
- b. 4 electrical tips and 16 pieces of LED lamps for treatment

- c. Daily skin care for body & face therapy. Stop time flowing through your skin
- d. Handy device is easy to carry around, treatment fulfill everywhere at anytime
- e. Enjoy skin care at home with in-office experience

### III Functions

The perfect synergy of any type of therapies

An innovative concept combining methotherapie and electroporation, we call methoporation it helps beauty essence into your skin easier 10 times or more than iontophoresis method. It is good for whitening skin, cellulite treatment. LED, RF, EMS, method also gives you many Anti-Aging benefit as skin rejuvenation, detoxification, lifting skin, lymphatic drainage, treatment of acne lesions, skin radition

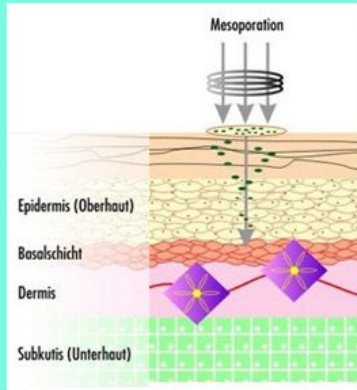


**“Mesoporation”** is a combination of mesotherapy and electroporation. In Europe it is widely used as mesotherapy without a needle. It achieves a high level of penetration into the skin. It is reported that Mesoporation’s penetration degree is about 10 times that of conventional therapy in a 1mm depth of skin and 27 times in 4mm. An amazing whitening and anti-aging effect is possible.




**For Treatment of Acne Lesions**


It works efficiently without damaging the skin by using an especially sensitive LED light developed exclusively for beauty equipment. It uses a blue LED with a wavelength of 415nm (±10nm).

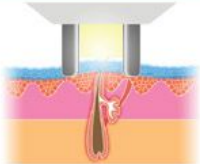

Contact is made with the skin and cooling begins.

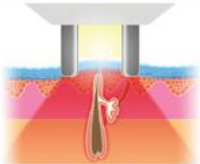

Optical energy pre-heats the target tissue.

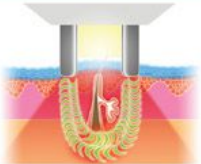

RF energy is drawn toward the heated target deep inside the dermis and both energies work to treat the target tissue.

### **a.Benefits For Face**

- (1). Lifting Skin
- (2). Whitening Skin
- (3). Treatment of Acne Lesions
- (4). Lymphatic Drainage
- (5). Skin Rejuvenation
- (6). Skin Relaxation

### **b.Benefits For Body**

- (1). Detoxification
- (2). Lymphatic Drainage
- (3). Skin Rejuvenation
- (4). Skin Relaxation

## **IV Specification**

a.Input Frequency:50-60Hz

b.Voltage:100-240V

c.Power

d.Consumption:10W

## **V Package Listing**

a.1 pc Main machine

b.1 pc Power Adapter

c.1 pc package box

## **VI How To Use**

### **a.Before use**

Please connecting the power cord with the machine,and clean your skin

### **b.Process**

(1).Please turn on the machine , chose the "Level 1 2 3 4." to chose different Level(1-4 From weak to strong),Don't turn it the strongest at first,You can adjust it from 1-4 to find the comfortable Level.

(2).Please put this device pressed to skin, but do not stay in the same place for too long,it need to move.( After ten minutes , the power will be shut down automatically)

(3).In the process , if you need to interrupt, stop, please press the LEVEL adjustment knob, then turn the power off .

## VII Operating Steps


Once 2 to 3 days, 15 to 20 times a course of treatment

### a.The leg(front):10 minutes

(1)From middle to two sides,it can active skin.



(2)From down to upper push, lift the skin.

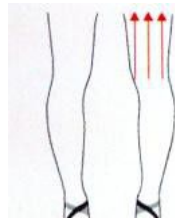


(3)From down to upper circling push.




### b.The leg(behind):10 minutes

(1)From down to upper push,lift the skin .

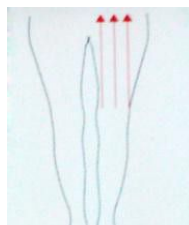




(2) From down to upper circling push.

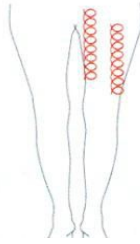


**c. The waist: 10 minutes**

(1) From down to upper push, lift the skin.

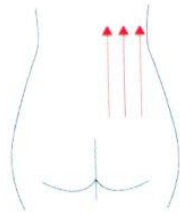

(2) From left to right push or from right to left .

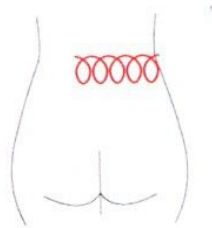

**d. The arm: 5 minutes**


(1) From down to upper push, lift the skin.



(2) From down to upper circling push, Stimulate the production of collagen.




(3) Outward from the arm looped, promote the effect of light, whitening, tender skin .

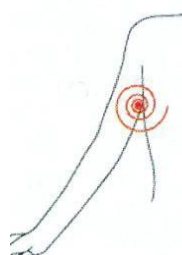

**e.The hip: 5 minutes**


(1) From down to upper circling push, active the skin.

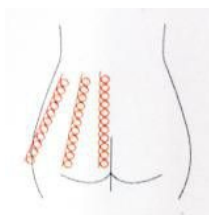

(2) From the crotch to the middle push,lift the hip, make Better hip shape .




(3)From down to upper push, lift the hip stimulate the production of collagen.

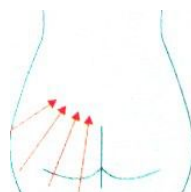


**f.The waist:10 minutes**


(1)push it as the line,active the skin on waist.

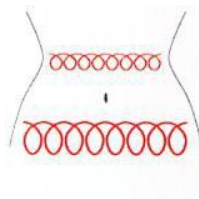

(2) In the middle of the waist, circle push around your belly button.

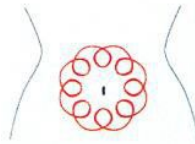

(3) From the waist of the lateral to the middle, Push it as pointed direction

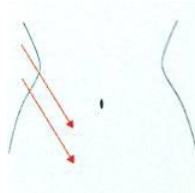

**g. The face: 10 minutes.**

Push it as the pointed direction.



**h. The eye: 5 minute    The lips: 5 minutes    The neck: 5 minutes**

Push it as the pointed direction.



## **VIII Forbidden Matter**

Please noticed that this item is not applied to people as following:

- a. People with heart disease or high blood pressure, or configured cardiac pacemaker
- b. Patients with acute inflammation, asthma, deep vein thrombosis, theyroncus, cancer
- c. Pregnant women
- d. People with hemorrhagic disease, trauma or who is bleeding
- e. Medical plastic parts, or parts with artificial in filling inside
- f. People with metal inside the body (include metallic tooth)
- g. Patient with skin inflammation or with edema
- h. People with an abnormal immune system
- i. Numb or insensitive to heat

## **IX Daily Inspection And Maintenance**

- a. If the head stick to a cosmetics or dirty, please use soft cloth gently wipe dry;
- b. Please don't tie line or winding, otherwise might well contact ;
- c. Thinner, volatile oil will be damaged, please don't use;
- d. Please do not placed in the sunlight, dust, more humid place (It is not waterproof);
- e. Please put it out of reach of children .

