




Operating Manual

Activity 845

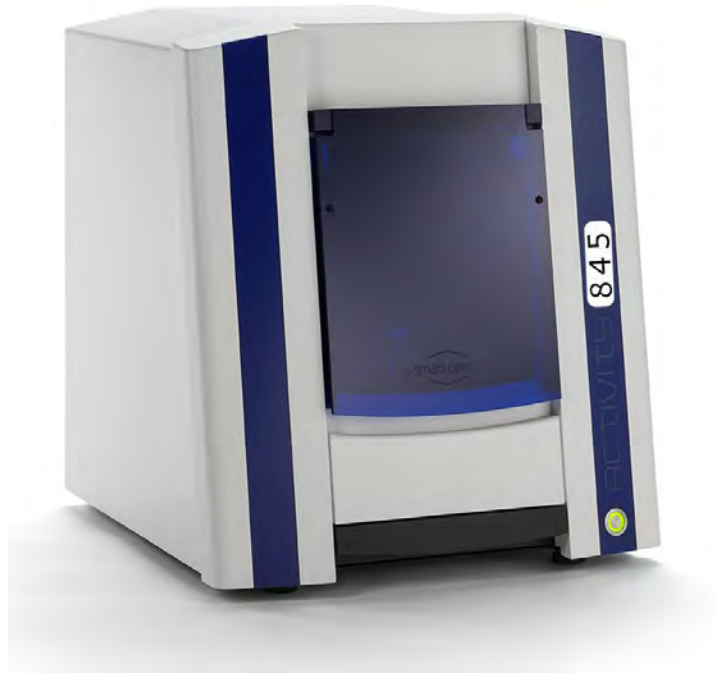A product of smart optics Sensortechnik GmbH





Activity 845

Table of contents

1. LIST OF FIGURES	4
2. ICONS	5
3. GENERAL SAFETY INFORMATION	6
4. TECHNICAL SPECIFICATIONS	6
5. CE DECLARATION OF CONFORMITY	8
6. GENERAL SAFETY INFORMATION	9
7. INSTALLATION	12
8.1 SELECTING THE INSTALLATION LOCATION	12
8.2 UNPACKING AND SCOPE OF DELIVERY	12
8.3 CARRYING POINTS	13
8.4 REMOVING THE SCANNER FROM THE PACKAGING	14
8. INSTALLATION	16
9.1 INSTALLING THE SCANNER	16
9.2 INSTALLING THE COMPUTER	17
9.3 INSTALLING THE SOFTWARE	17
9.4 IMPORTING CALIBRATION DATA	23
9. CALIBRATION PROCEDURE	26
10.1 AXIS CALIBRATION	26
10.2 3D CALIBRATION	28
10. BASIC INFORMATION ABOUT THE DEVICE	30
11.1 FUNCTIONING OF THE ACTIVITY845	30
11.2 INTERIOR ACTIVITY845	34
11.3 OBJECT HOLDER	35



<u>11. SCANNING PROCEDURE</u>	38
12.1 CLAMPING THE MODEL	38
12.2 CREATING A PRESCAN	38
12.3 CASE STUDY: BRIDGE WITH BITE REGISTRATION	41
12.4 CASE STUDY: COMPLETE JAW SCAN	50
12.5 CASE STUDY: OCCLUSION SCAN	54
<u>12. SYMBOLS</u>	62
13.1 THE SYMBOLS	62
13.2 THE 3D VIEWER SYMBOLS (OBJECT VIEW)	65
13.3 ACTIVITY MENU OPTIONS	66
13.4 OPTIONS → SETTINGS → MATCHING	67
13.5 OPTIONS → SETTINGS → GENERAL	68
13.6 OPTIONS → SETTINGS → INSTALLATION	72
<u>13. MAINTENANCE AND SERVICING</u>	73
<u>14. FAULTS AND REPAIRS</u>	73
<u>15. ENVIRONMENTAL FACTORS AND DISPOSAL</u>	73
16.1 PACKAGING	73
16.2 DISPOSAL	74
<u>16. IMPRINT</u>	74

1. List of figures

FIGURE 1 SCANNERSIDE VIEW	15
FIGURE 2 SCANNER HOLDING POINT	15
FIGURE 3 : SCANNER HOLDING POINT	15
FIGURE 4 MAIN SWITCH ON THE FRONT OF THE DEVICE	16
FIGURE 5 SETUP/INSTALLATION.....	18
FIGURE 6 LANGUAGE SELECTION	18
FIGURE 7 SETUP WIZARD	19
FIGURE 8 SELECT INSTALLATION LOCATION.....	20
FIGURE 9 START MENU.....	20
FIGURE 10 CREATING A DESKTOP SYMBOL.....	21
FIGURE 11 INSTALLATION.....	22
FIGURE 12 INSTALLATION.....	22
FIGURE 13 COMPLETING THE INSTALLATION	23
FIGURE 14: NO CALIBRATION DATA.....	24
FIGURE 15: FOLDER SEARCH	24
FIGURE 16: SELECTED FOLDER.....	25
FIGURE 17: IMPLEMENTING THE SENSOR DATA	25
FIGURE 18: CALIBRATION MODEL WITH SPACER PLATE IN ACTIVITY845.....	26
FIGURE 19: SCAN SOFTWARE CALIBRATION	27
FIGURE 20: SCAN SOFTWARE CALIBRATION	27
FIGURE 21 INDIVIDUAL CALIBRATION BLOCK	28
FIGURE 22: CALIBRATION MODEL REGISTRATION.....	29
FIGURE 23: COMPLETING THE CALIBRATION	29
FIGURE 24: DEVICE FRONT A845	32
FIGURE 25: DATA PLATE.....	33
FIGURE 26: ROTATING-SWIVELLING UNIT	34
FIGURE 27: ROCKER/ROTATING-SWIVELLING UNIT ACTIVITY845	34
FIGURE 28: INDIVIDUAL PARTS OF ACTIVITY 845	36
FIGURE 29: OBJECT HOLDER SYSTEM	37
FIGURE 30: OBJECT HOLDER SYSTEM WITH MODEL.....	37
FIGURE 31: "NEW PROJECT"	38
FIGURE 32: "NEW MEASUREMENT".....	39
FIGURE 33: 2D SCAN AND ENTRY SCREEN	40
FIGURE 34: PROMPT NOTIFICATION	42
FIGURE 35: DATA SET BEFORE "MATCHING"	43
FIGURE 36: "RESCAN" MODE.....	45
FIGURE 37: "RESCAN".....	45
FIGURE 38: SAVE PROJECT.....	46
FIGURE 39: "MATCHING" PROCESS.....	47
FIGURE 40: 3D STL DATA SET	47
FIGURE 41: "BITE REGISTRATION" SYMBOL.....	48
FIGURE 42: INFORMATION WINDOW	48
FIGURE 43: SQUEEZE BITE ALIGNMENT PRESCAN	49
FIGURE 44: CUT BITE REGISTRATION.....	49
FIGURE 45: CUT BITE REGISTRATION.....	50

FIGURE 46: COMPLETE JAW 2D SCAN	51
FIGURE 47: COMPLETE JAW ALIGNMENT PRESCAN	52
FIGURE 48: STARTING THE SCAN PROCESS	52
FIGURE 49: COMPLETE JAW 3D SCAN	53
FIGURE 50: 3D-STL DATA SET	54
FIGURE 51: TABS	54
FIGURE 52: FIXED JAW.....	56
FIGURE 53: ADD DENTAL SCAN.....	56
FIGURE 54: NEW MEASUREMENT	57
FIGURE 55: ANTAGONIST ALIGNMENT	57
FIGURE 56: PRESCAN	58
FIGURE 57: VESTIBULAR SCAN.....	58
FIGURE 58: "ARTICULATION" ICON	59
FIGURE 59: SELECT LOWER JAW	59
FIGURE 60: SELECT UPPER JAW.....	59
FIGURE 61: LOWER JAW ALIGNMENT.....	60
FIGURE 62: "NEXT"	60
FIGURE 63: UPPER JAW ALIGNMENT	61
FIGURE 64: ALIGNED UPPER AND LOWER JAWS	61
FIGURE 65: ARTICULATED STL DATA SET	62
FIGURE 66: ACTIVITY MENU OPTIONS	66
FIGURE 67: SETTINGS - MATCHING	67
FIGURE 68: SETTINGS - GENERAL.....	68
FIGURE 69: SETTINGS - INSTALLATION.....	72

2. Icons

	This symbol indicates warning information
	Useful hints are marked with a light bulb in the operating manual

3. General safety information

Proper use:

The Activity845 scanner is intended for use in performing optical, three-dimensional measurements of human jaw models.

It is possible to measure articulated models in order to simulate masticatory movements using corresponding software, as in an articulator.

Anything other than the proper use is strongly advised against, as improper use may cause damage or injury.

4. Technical specifications

	Activity845
Axis	2
Dimensions	440 mm x 480 mm x 430 mm (WxHxD)
Weight	32,2 kg
Connection voltage	100-240 V AC, 50/60 Hz
Power consumption	80 W
Protection class	IP11
Permissible temperature range	18-30 °C
Measureable objects	Plaster models of teeth
Minimum clamping range of the removable object holder	40 mm
Maximum clamping range of the removable object holder	70 mm
Height adjustment	//

Measuring time for the measurement of	Single stump: < 55 sec./unit 3- bridge: < 90 sec. Complete jaw: < 90 sec.
Resolution	Basic resolution of the 3D sensing head:62,5 µm 80 x 60 x 85 mm (x,y,z)
Accuracy	+/-10 µm standard deviation measured at test specimen, determined via min. 50 measuring points with 10 repetitions each
Output data format	STL
Interfaces	USB
Scope of delivery	Scanner, power cord, 2x USB cable, CD with operating software and calibration data, operating manual, calibration block
System requirements	Operating system Windows 7 Ultimate 64® Bit, Intel core 2 Quad CPU Q 9550 2.83 Ghz or higher, RAM 8 GB, high-performance 3D graphics card with min. 1GB RAM e.g. GeForce, 320 GB hard disk



Sensortechnik GmbH
Lise-Meitner-Allee 10
44801 Bochum, Germany

Declaration of CE conformity

according to EU-regulation 2006/42/EG

We declare that the device identified below complies with the requirements of the EU guideline which regard to safety and physical health requirements both in concept and construction put in circulation.

This declaration becomes invalid in case of an unauthorized change of the device.

Device description: Optical 3D scanner

Device type: **Activity 845**

EU guidelines applicable: machine guideline (2006/42/EC)
 low voltage guideline (2006/95/EEC)
 EMC guideline (2004/108/EEC)

Harmonized standards applied:

EN 1050, EN 12100-1, EN 12100-2, EN 61000-6-1, EN 61000-6-3

The CE label was used first for this product in 2013.

Document prepared by: Jörg Friemel

Signature:

A handwritten signature in black ink that reads "Jörg Friemel". The signature is written in a cursive style with a large 'J' and 'F'.

Bochum, 12.05.2013

smart optics
Sensortechnik GmbH
Lise-Meitner-Allee 10
44801 Bochum
Tel.: 0234/29 82 8 - 0 Fax: -20

6. General safety information

When setting up, commissioning and operating the scanner, you must observe the following safety information at all times.

1.

The device may be operated only by adequately trained staff who have knowledge concerning the proper use of the device as well as the safety information cited here.

2.

The device is intended exclusively for use in dry, enclosed rooms.

3.

The device may be operated only on a stable base (tabletop, bench, etc.) that has a load bearing capacity which is at least twice the net weight of the scanner as stated in the technical specifications.

The base must be fixed to a building wall or have a braced support frame that is suitable for absorbing vibrations occurring during operation of the scanner.

Besides the pure load-bearing capacity, the sturdiness and stability of the base is crucial for safe operation of the scanner.

4.

The information on transport and storage must be observed when lifting and transporting the scanner.

5.

Please do not place any objects on the scanner, as there is a risk that these can fall down owing to the surface structure as well as the resultant vibrations during scanning operation.

6.

The “**Activity 845**” scanner has been developed and manufactured in accordance with the applicable safety standards and with the greatest possible care. In spite of this, the risk of an electric shock, overheating or fire due to technical defects in individual components cannot be totally excluded.

You should therefore switch the device off when not using it for longer periods and when it is left unattended. This also benefits the environment, as it means less electrical energy is consumed.

7.

The scanner contains a motion unit comprising three electrically powered axes of rotation that serve to position the object during scanning.

To enable the object to be clamped in a definite position as well as a faster workflow during scanning, the movement axes are also kept in position electrically when the flap is open.

An electronic protection mechanism prevents the motors from starting up unintentionally or continuing to operate when the flap is opened.

As there is a theoretical possibility of the protection mechanism failing, this entails a residual risk which you can counter by complying with the following **rules of use**:

7.1.

Do not open the front flap during a scanning procedure until after a corresponding prompt or the end of a scanning procedure has been displayed on the screen.

7.2.

If you determine that the scanner is acting unusually, such as an uncontrolled or continuous rotation of one or more movement axes, switch off the scanner at the main switch before opening the flap.

If this unusual behaviour continues after switching on again and restarting the software, do not use the scanner any longer. The device must be marked as defective and prevented from being put into operation again until the fault is rectified.

7.3.

Do not reach inside the scanner if one or more movement axes move despite the flap being open. **Switch the scanner off immediately.**

The scanner may no longer be used in this case and must be marked as defective and prevented from being put into operation again until the fault is rectified.

7.4.

During scanning, never reach below the area between the rocker and case that is marked by warning signs.

This also applies if the rocker is swivelled upwards to the rear.

Before removing any objects that have fallen into the scanner or before cleaning the scanner compartment, you must switch off the device at the main switch.

8.

If you notice that the scanner is damaged or defective, you must label it as defective and make sure that it is not used until it has been repaired.

7. Installation

8.1 Selecting the installation location

Before installing the scanner, you should select a suitable installation location.

A suitable, stable base (bench, tabletop etc.) should be chosen as an installation location.

If you install the computer required for operating the scanner underneath the worktop, the available worktop area should be at least **1.1 m x 0.75 m (front width x depth)**. If you plan to also place the computer on the worktable, the area must be enlarged corresponding to the dimensions of the computer.

The selected workplace should, if possible, face away from windows or strong artificial light, as excessively strong external light can lead to undesirable reflections on the monitor that can interfere with the functioning and results of the optical scanner in extreme cases.

8.2 Unpacking and scope of delivery

Inspect the external packaging for visible damage as soon as you take delivery of the scanner. If you notice that the packaging is damaged, inform the delivering forwarder and your specialist dealer immediately.

The device is delivered in a strong outer box on a wooden pallet. After removing the straps, open the upper folding lid of the box. The box contains a foam piece that protects the scanner, which is packed below during transport.


First, pull the foam piece upwards out of the box. You can then lift off the whole box upwards.



The net weight of the scanner is approx. 32,2 kg.

Always observe the following instructions to lift the scanner properly out of the packaging and position it at its location of use.

- Two people must lift the device out of the packaging and place it on the prepared, intended workplace.
- Be sure to transport the scanner (while still on the pallet) as closely as possible to the workplace to avoid longer transport distances.

	<p>The scanner accessories can be found at the base of the packaging. Check that the scope of delivery is complete.</p>
---	---


Scope of delivery	Activity845
1 scanner	X
2 USB cables	X
1 hex key	X
1 user manual in PDF format	X
1 calibration object	X
1 software installation CD including calibration data	X
1 removable object holder	X
1 measuring range template	
1 spacer plates 10mm	

8.3 Carrying points

Carrying points are provided to move the device. Please note these specifications.

Stand so that one person each is positioned to the left and right of the scanner. Take hold of the scanner with one hand at the bottom front of the scanner and the other hand at the top rear. In this way, the scanner can be safely tilted slightly to the rear.

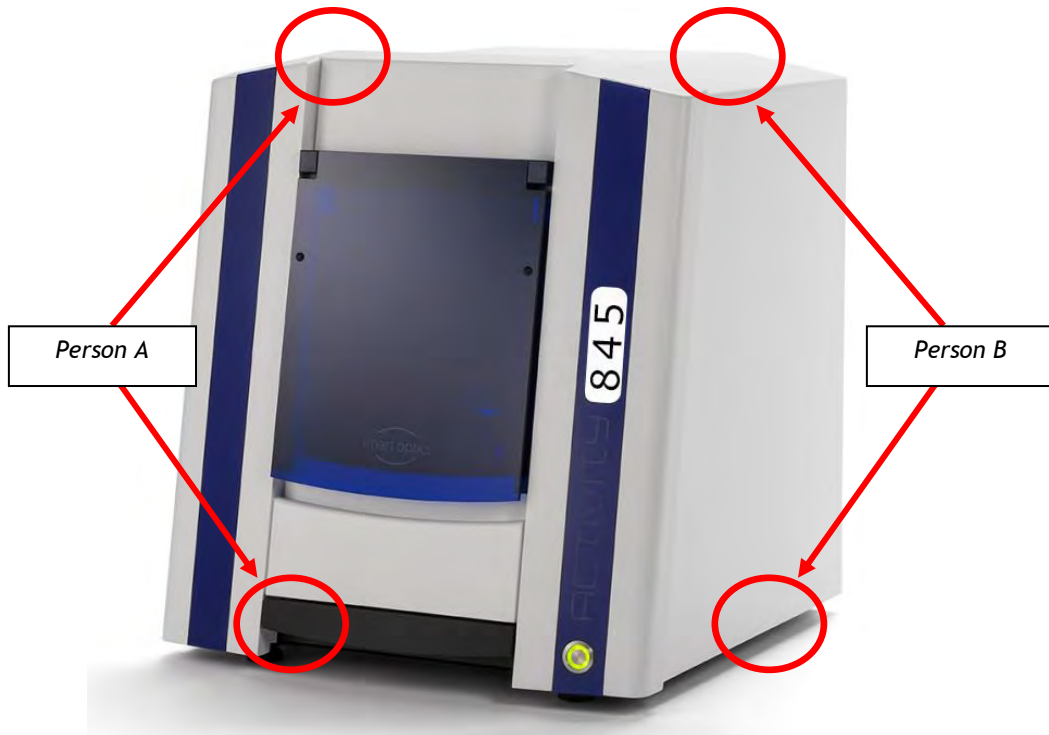Carrying points, front



Never carry the scanner with the flap open, and do not carry the scanner using the flap itself.

8.4 Removing the scanner from the packaging

1. To remove the scanner from the packaging, one person stands to the left of the scanner and the other person stands to the right.
2. Each person must take hold of the scanner with one hand on the upper support point. Next, tilt the scanner slightly to the rear until you can grip underneath the scanner at the front with your other hand. (Figures 2 and 3)
3. Working together and at the same time, lift the scanner out of the packaging and place it at the installation site.

4. When doing so, make sure that you have free access to the workplace and do not trip over the packaging.



Remove the upper part of the packaging



Figure 1 scannerside view



Grip and tilt the scanner at the upper holding point (one holding point each side for each person).

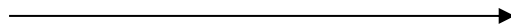


Figure 2 Scanner holding point



Lift the scanner at the same time and place it at the installation site.



Figure 3 : Scanner holding point

8. Installation


9.1 Installing the scanner

Please make sure that the mains switch is at the "0" position before installation. (The position of the main switch maybe different)



Figure 4 Main switch on the front of the device

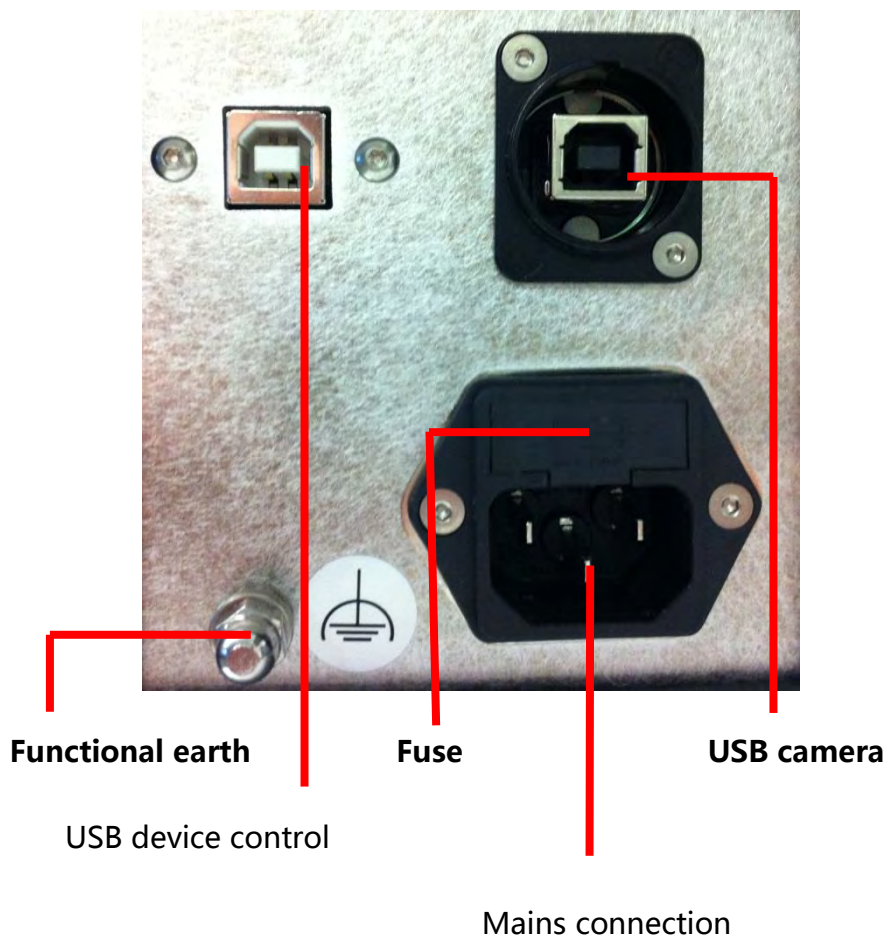1. Connect the device to the USB cables using the **USB ports** on the rear of the device.



2. Connect the other ends of the USB cables to USB connections on the **rear** of your computer.
3. Connect the scanner to a power source via the mains connection on the rear.
4. Now switch the scanner on at the mains switch.
5. Now first switch the scanner off again and continue by installing the computer and operating software.

9.2 Installing the computer

Install the computer following the installation instructions of the computer manufacturer.

9.3 Installing the software



As each scanner has been calibrated individually at the factory, the data on the installation CD is valid only for the corresponding scanner. This calibration data contains a code that can be enabled only by the corresponding scanner. Therefore, compare the serial number indicated on the data plate of the scanner with the serial number indicated on the CD to make sure they match **before** installation.

If you have mistakenly installed an incorrect version of the software, you must uninstall the incorrect version first.

If you are installing the Activity software on your computer for the first time, use the installation CD provided with the scanner.

A window will open after inserting the CD. The "Setup" icon opens and must be selected to start the installation.

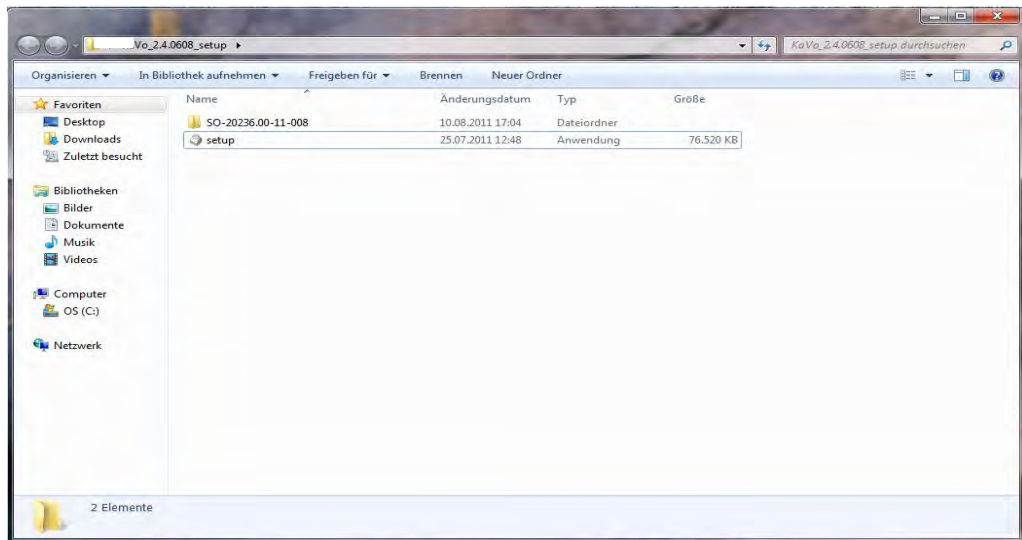


Figure 5 Setup/Installation

Now select the required language (Fig. 6) and confirm your selection with the "Next" icon (Fig. 7).

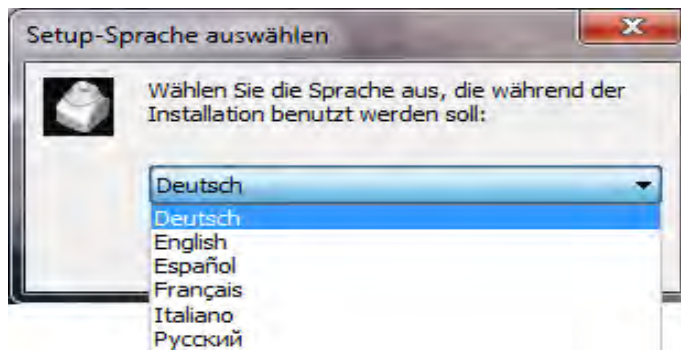


Figure 6 Language selection



Figure 7 Setup Wizard

Define the installation location for the software by using the "Browse" button to define the path. The installation program suggests a standardized path that you can use. Having selected the path, please confirm this with "Next".

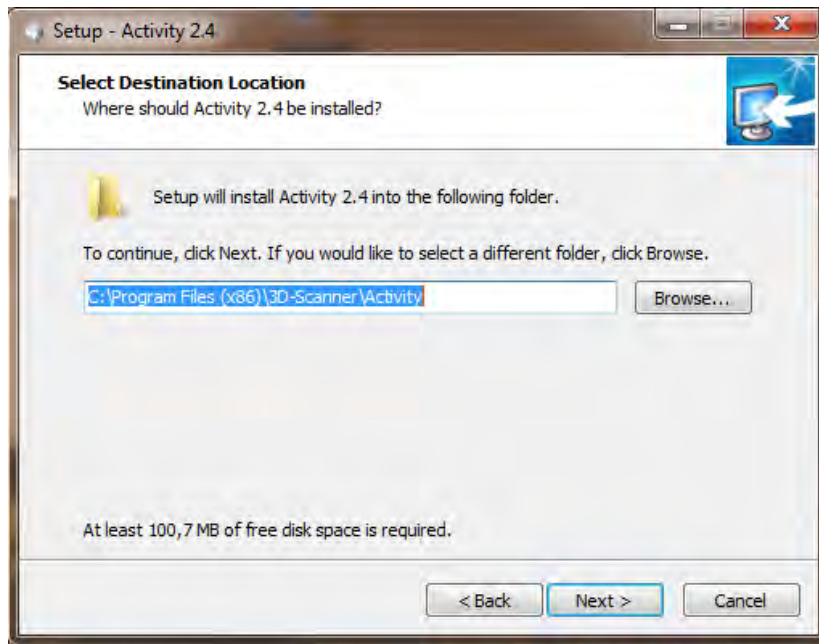


Figure 8 Select installation location

Define the software name in the "Select Destination Location Folder" window. Then click on "Next".

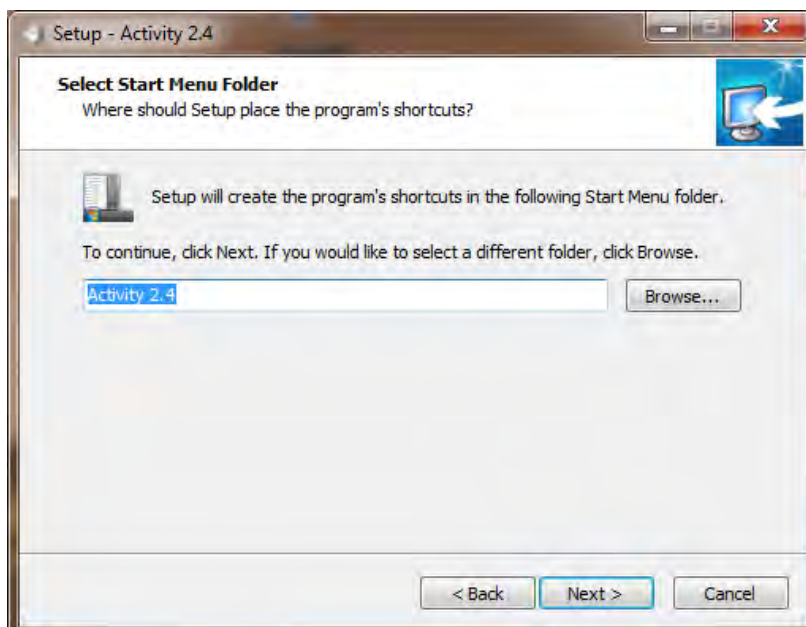


Figure 9 Start Menu

Select whether you require a Desktop symbol and confirm with "Next".

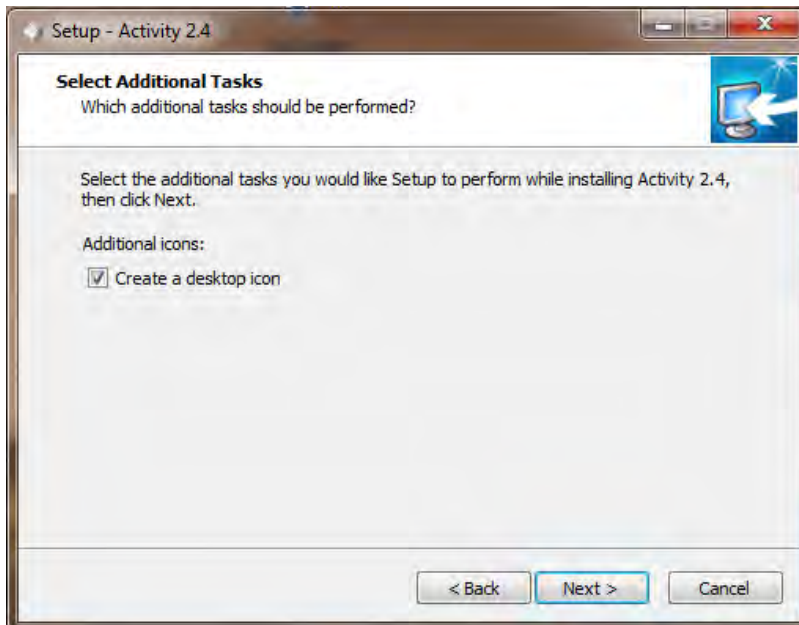


Figure 10 Creating a Desktop symbol

Before starting the installation, check your selected settings and chose "Next". The installation starts immediately (Fig. 11).

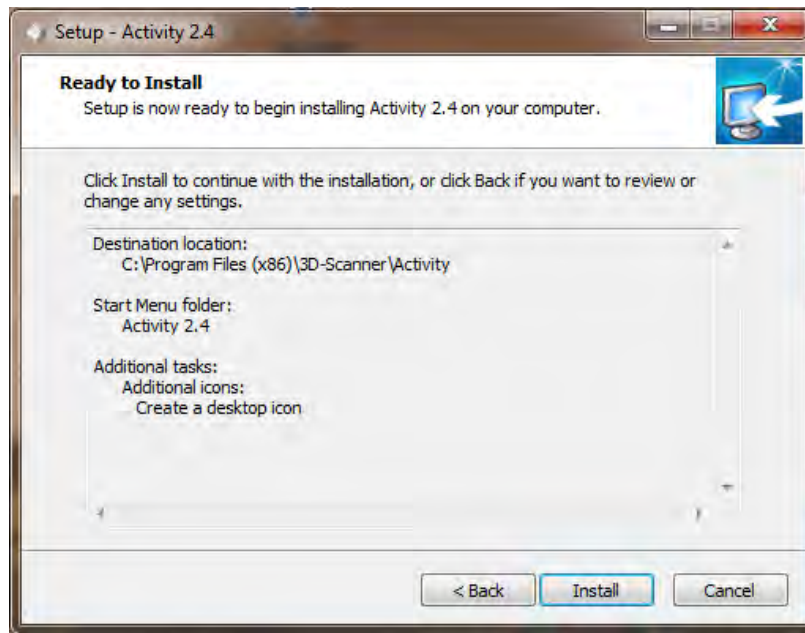


Figure 11 Installation

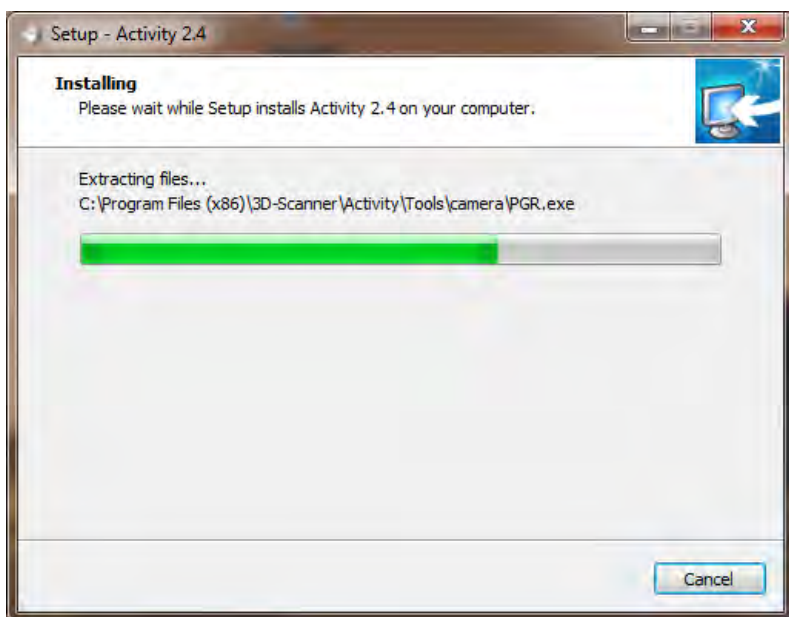


Figure 12 Installation

Confirm the installation of the driver with "Next". Then complete the installation with "Finish". The drivers are now installed and the process is completed.



Figure 13 Completing the installation

9.4 Importing calibration data

The scanner-specific calibration data must be imported during the initial installation.

Start the Activity software via the Desktop icon or from the Start Menu.

This is followed by a software message (Calibration directory not found). Click "OK" to start the "Installer Tool".

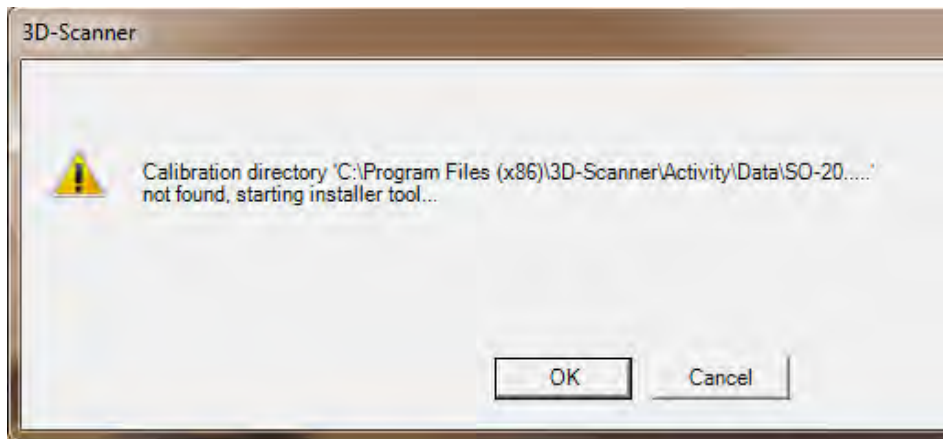


Figure 14: No calibration data

This tool makes it easier to install the required sensor data.

Click "Search" to select the file path (**SO-202....**).

When doing so, select only the parent folder and confirm with "OK".

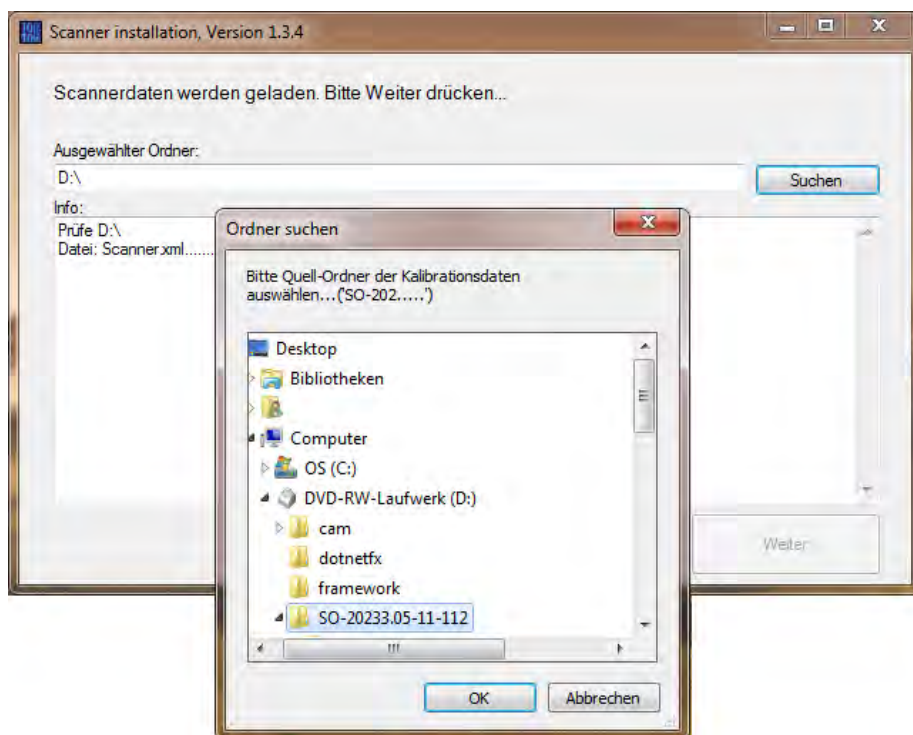



Figure 15: Folder search

Click "Next".

The calibration data is imported into the Activity software.

9. Calibration procedure

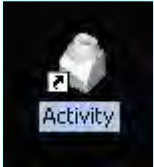10.1 Axis calibration




Start the scan software by clicking on the Activity icon located on your Desktop.

Insert the **spacer plate** in the Activity850 scanner. Clamp the calibration model on the removable object holder and place it in the scanner (Fig. 18).



The Activity845 can be calibrated without a spacer plate.

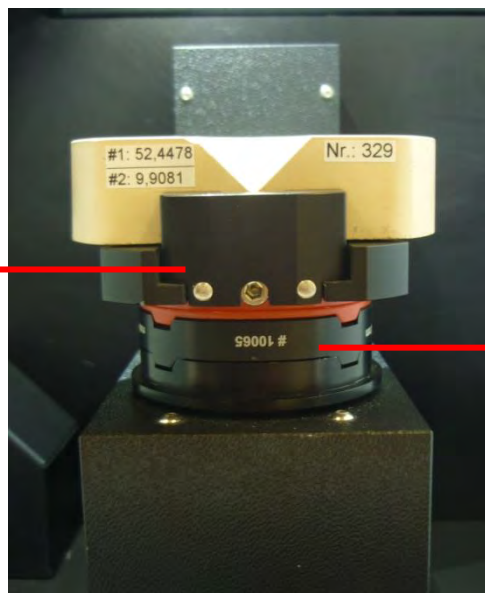Object base



Spacer plate

Figure 18: Calibration model with spacer plate in Activity845

The path **Options**→ **Service**→ **Axes Calibration** contains the start command for calibrating the axes.

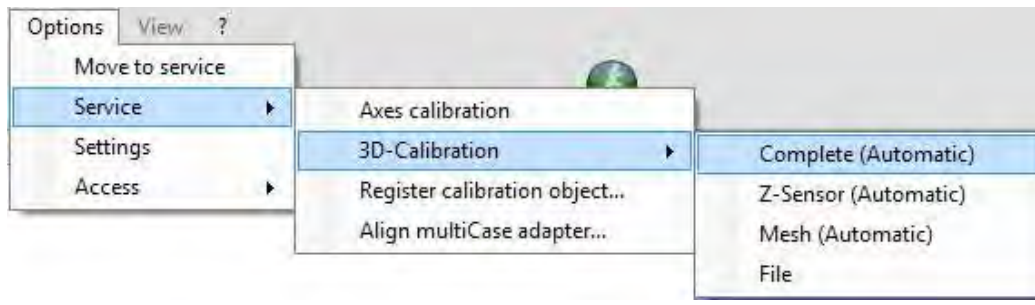


Figure 19: Scan software calibration

The scanner now performs an automatic axis calibration. The following message will appear after the axis calibration is complete:

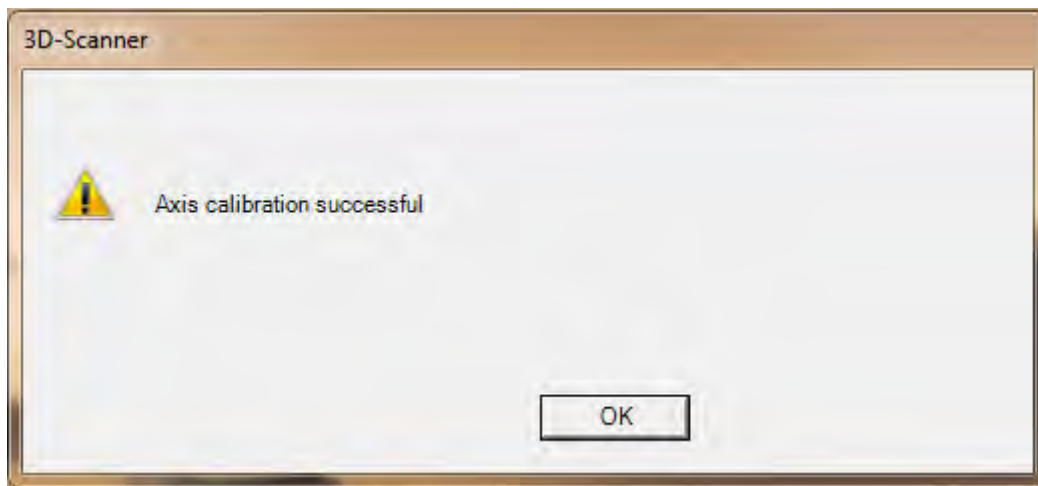



Figure 20: Scan software calibration

	<p>If the axis calibration fails, please check the height alignment of the model in the measuring field. (See Figure 18)</p>
---	--



Each calibrating block is measured industrially.
These individual values can be found on the calibration model.

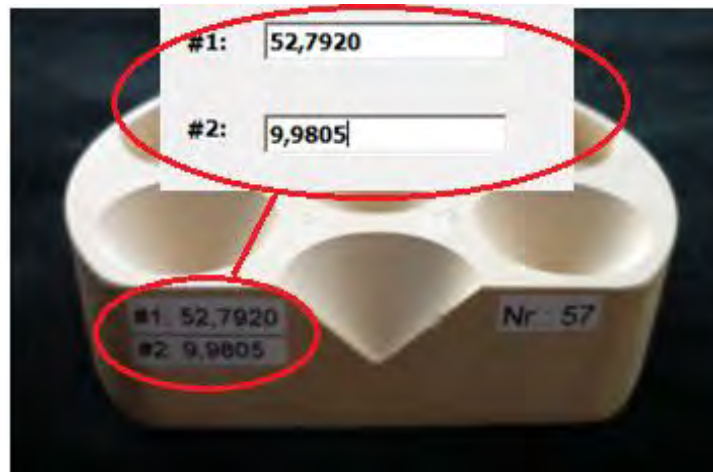


Figure 21 Individual calibration block

These values must be saved in the software after reinstalling the software or replacing the calibration model.

Start this procedure as follows:

A window opens under **Options** → **Service** → **Calibration object registration** (Fig. 22).

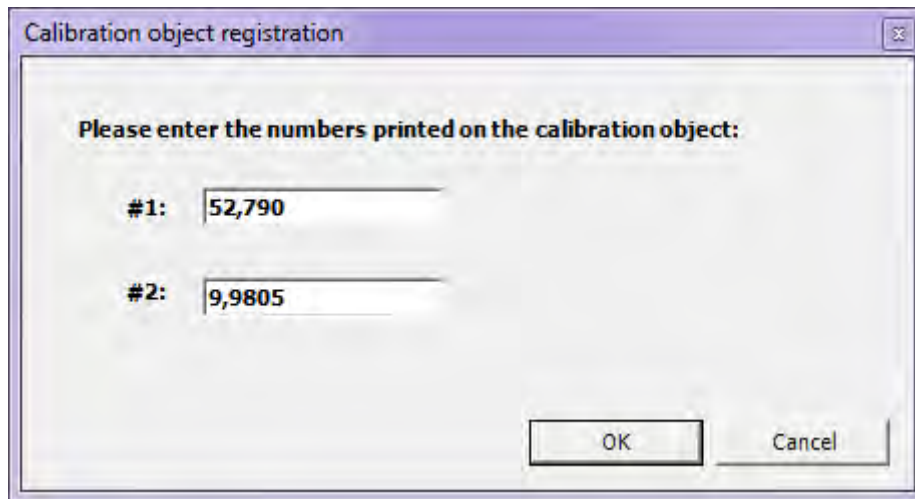


Figure 22: Calibration model registration

Individual values are indicated on the rear of the calibration model. Please enter these in fields #1 and #2.

Confirm the procedure with "OK".

You can start the 3D calibration only by selecting the start command under the item **Options→ Service→ 3D calibration → Complete**.

The system now performs a 3D calibration; this can take several minutes.

The following message will appear at the end of the calibration (Fig. 23).

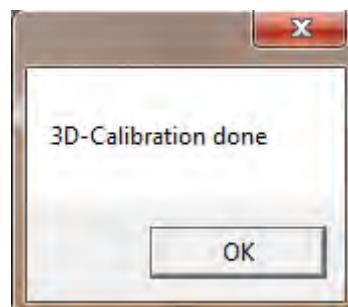



Figure 23: Completing the calibration

The 3D calibration was successful. Confirm with "OK".

	<p>Repeat the calibration about every four weeks and each time the scanner is transported in order to ensure the accuracy of the scan results.</p>
---	--



If the **axis calibration** fails, please check the height alignment of the model in the measuring field. (See page 16, Fig. 18)

10. Basic information about the device

11.1 Functioning of the Activity845

The optical scanner with the designation "**Activity845**" is used for the three-dimensional measuring of jaw models, in orthodontic and prosthetic applications.

The most important components of the scanner are the **3D sensor** and the **positioning unit** with object holder.

The **3D sensor** consists of a **camera** and a **projector**. The projector throws a light pattern onto the object to be measured at an oblique angle of incidence; the object is simultaneously monitored with the camera.

Since the camera cannot capture the entire object if the object is complex in shape, as in the case of a jaw model, the object is mapped with a series of individual scans from various perspectives.

The software then combines the individual scans to create a complete data set.

The positioning unit positions the object to be measured in the required perspectives opposite the 3D sensor.

It contains motorized axes for this, which can rotate and tilt the model on the base relative to the 3D sensor.

The entire scanning procedure is extensively automated so that you only have to carry out the essential operating steps.

A 3D scan is created as follows:

The model to be measured is mounted on the object holder:

The model is fixed on the object holder for this.

Prescan:

The prescan compiles a two-dimensional image of the clamped model.

This image serves as a basis for the following scan definition.

Scan definition:

The scan definition defines the positions in which the model is to be scanned.

The scan definitions are formulated using the prescan image after the user highlights the tooth positions on the screen.

Individual measurements:

After completing the scan definition, the software automatically calculates a scan program for recording the positions specified by the operator and starts these. The positions provided in the scan program are started up by the positioning unit and a 3D measurement is performed.

When scanning adjacent tooth stumps, it is usually necessary to measure the stumps in isolation without the adjacent stumps, as otherwise the interdental area cannot be captured by the 3D sensor.

The software will therefore prompt you to present the relevant tooth stumps individually and freely, i.e. to remove all other stumps from the saw cut model.

Matching and data storage:

The recorded individual measurements are grouped together in the "Matching" mode.

The finished data set is then saved on the hard disk of the computer.

Visualization:

When the matching process is complete, the scan result is displayed on the monitor so that you can directly assess the result visually.

The most important components and operating controls of the scanner are explained below.

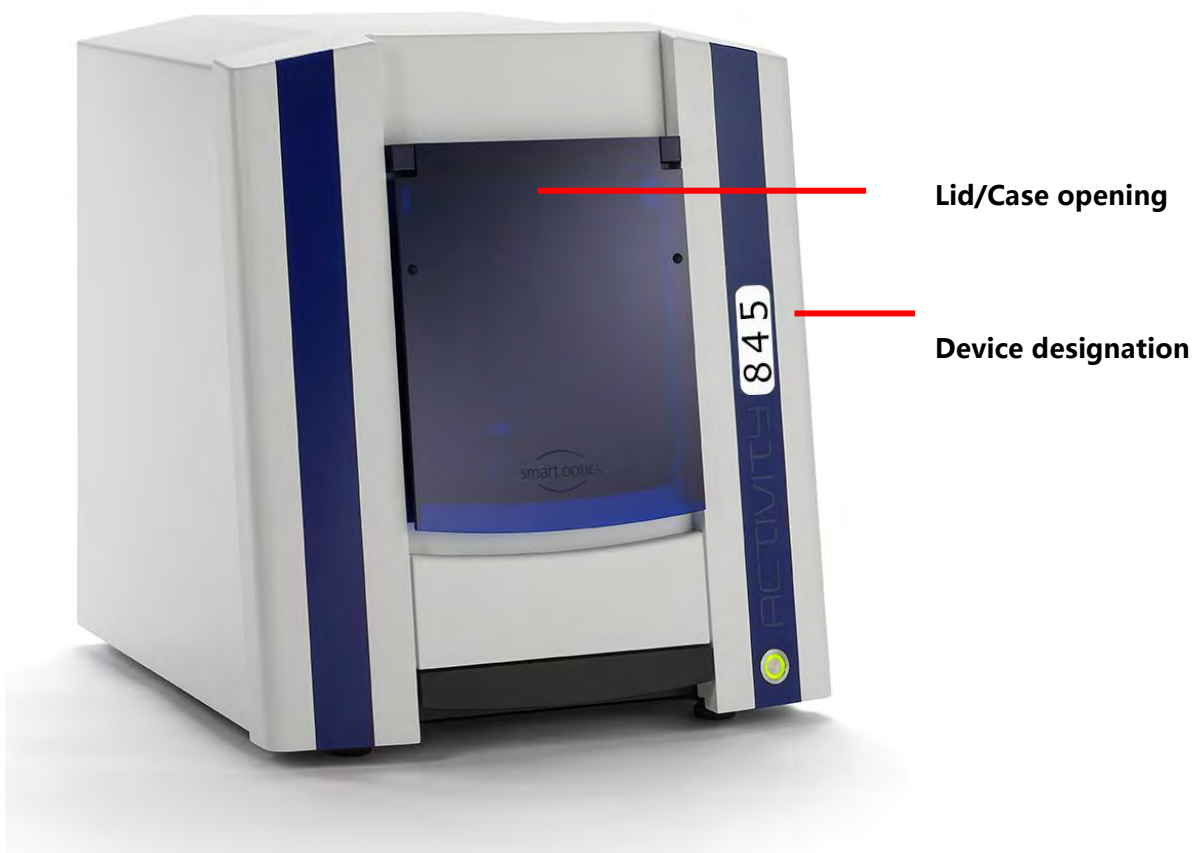


Figure 24: Device front A845

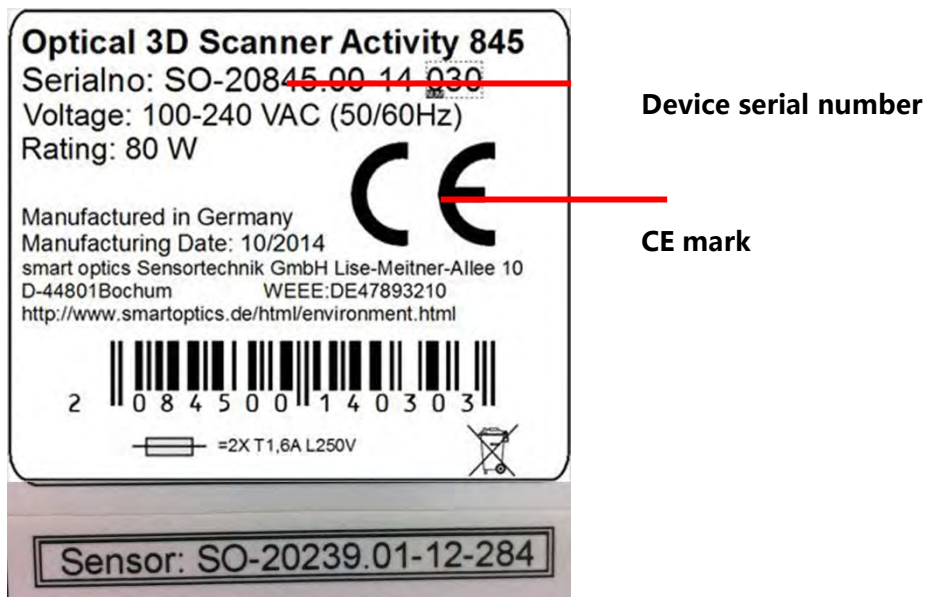


Figure 25: Data plate



Use only the USB cables included in the scope of delivery, if possible. These have undergone multiple testing at our company in conjunction with your scanner. Communication problems can result between the scanner and PC if cable lengths over 2 metres are used.

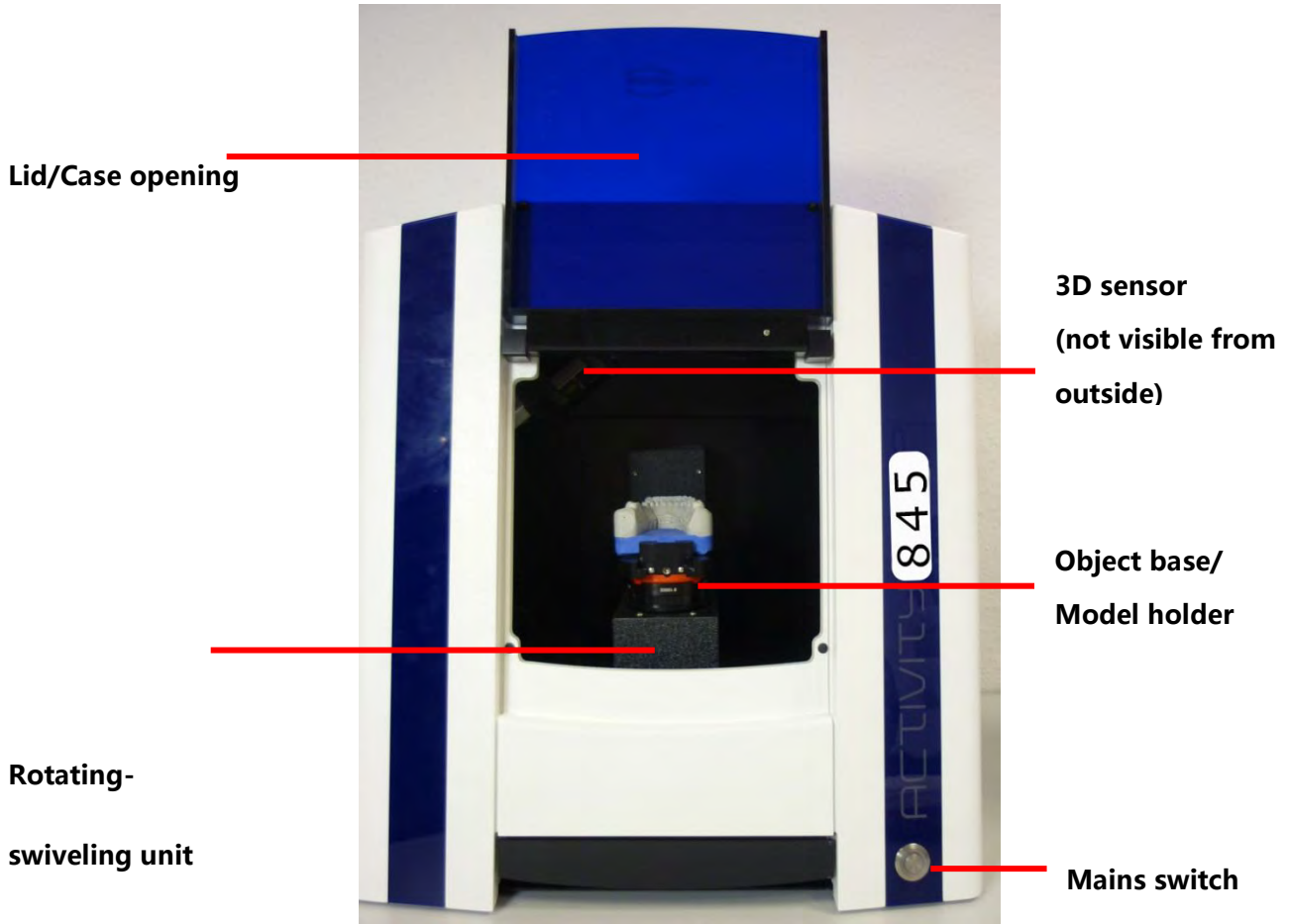


Figure 26: Rotating-swivelling unit

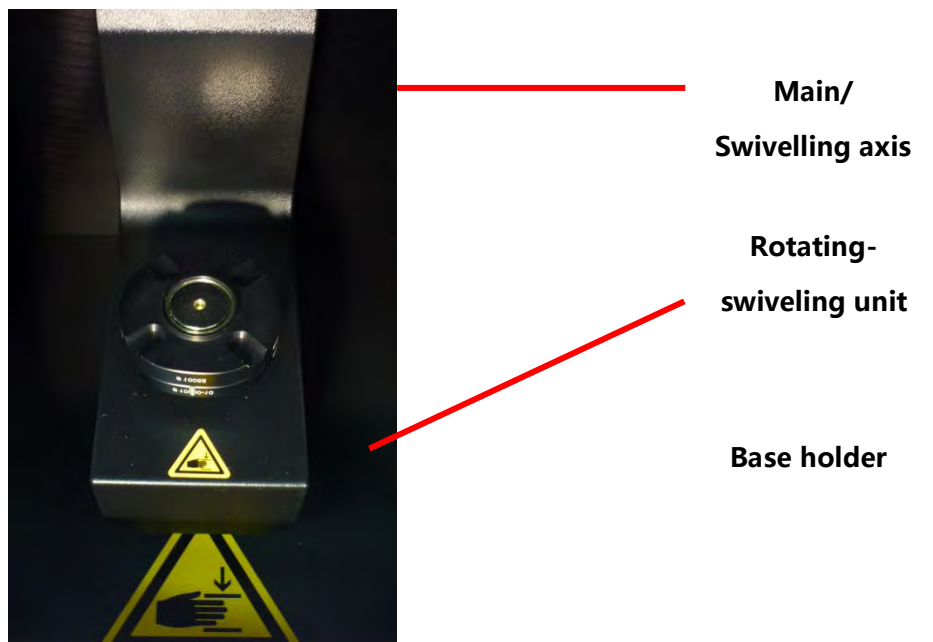


Figure 27: Rocker/Rotating-swivelling unit Activity845

Positioning unit (inside the scanner):

The positioning unit contains the fixed base holder on which the object to be measured is fixed as well as a rotating and swivelling unit powered by an electric motor.

Rotating and swivelling unit/Object holder:

The rotating unit allows you to position the object to be measured in any position relevant to the camera of the 3D sensor located on the object holder via the electric motor.

The swivelling unit allows you to swivel the entire setup comprising rotating unit and object holder and hence record side views of the object to be measured.

Interior lighting:

The interior lighting switches on automatically when the flap is opened. If the flap remains open for longer than 5 minutes, the lighting will switch off automatically.

The interior lighting switches off automatically when the flap is closed.

11.3 Object holder

Individual parts of the object holder

The scope of delivery for the object holder of the **Activity845** consists of the following components:

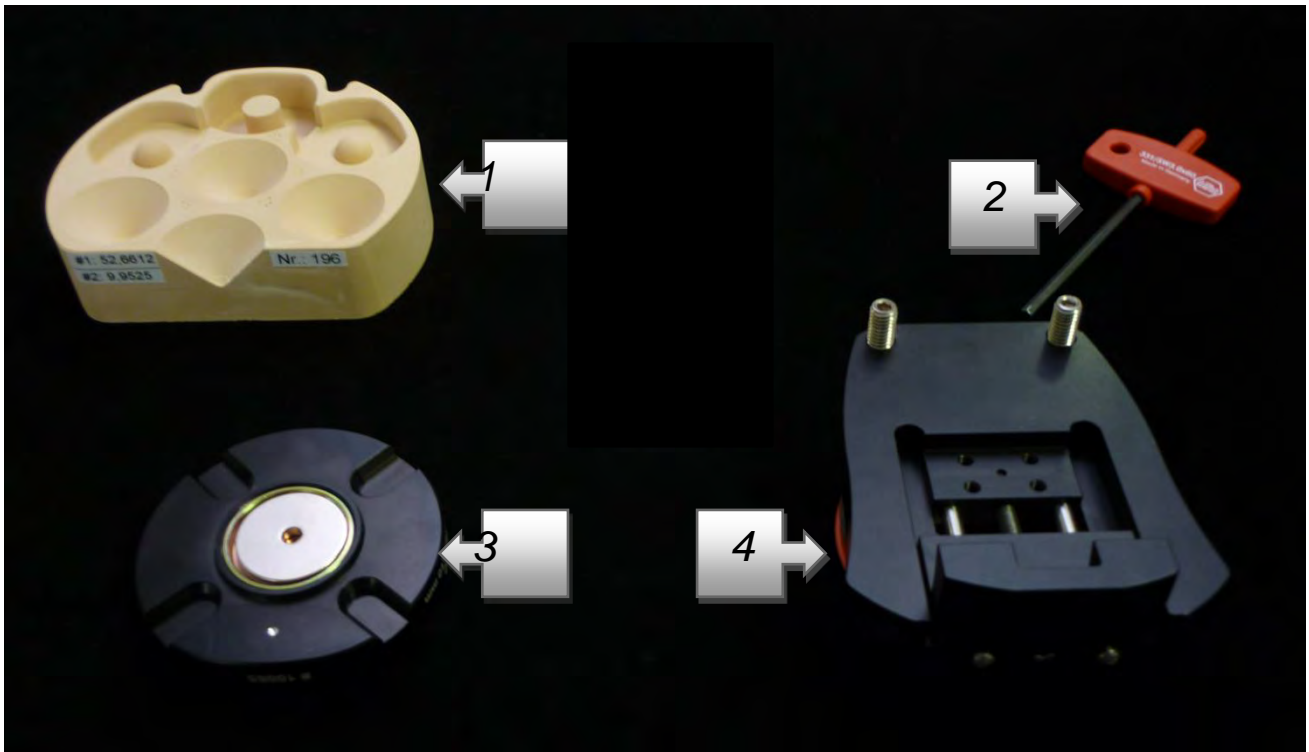


Figure 28: Individual parts of Activity 845

Key to the scope of supply of the object holder system

1. **Calibration object** for axis and 3D calibration.
2. **Hex key** for clamping or loosening the plaster model on the model holder.
3. **An additional spacer plate with magnet** for height alignment.
4. **Removable object holder** for non-articulated jaw models. The model is fastened on the object holder using a clamp. This clamp is fastened or released using the hex key (2).

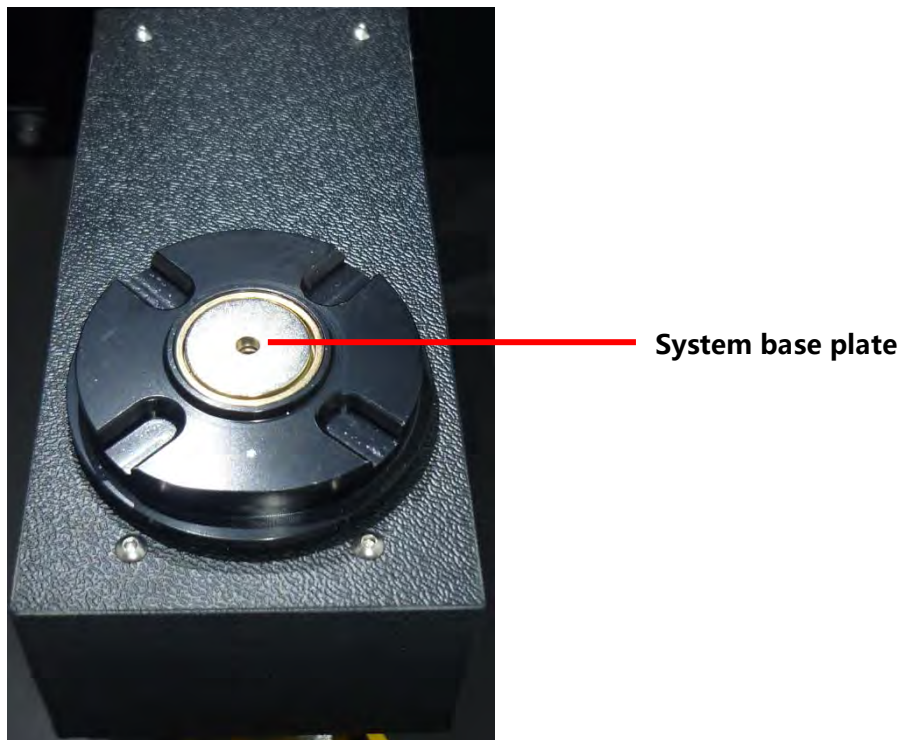


Figure 29: Object holder system

In order to ensure an optimum alignment of the model to be scanned in the measuring range, a maximum of one additional spacer plate can be mounted on this. The object holder can now be fixed on the base or spacer plate.

Example: Standard model holder without spacer plate



Figure 30: Object holder system with model

11. Scanning procedure

12.1 Clamping the model

Clamp the model in the removable object holder and align the height using the spacer plate (see Figure 30).

12.2 Creating a prescan



Depending on the scanner and CAD software combination, it is not necessary to fill out the project information thanks to a special interface, as this is automatically undertaken by the CAD software (e.g. with dentCreate!®). The scanner uses software messages to guide you step by step through the entire scanning procedure.

Switch the scanner on at the mains switch. The scan software opens after starting by double clicking on the "**Activity**" icon on your Desktop. An automatic referencing of the axes takes place. The project entry screen appears:

New Project

Laboratory: Include in name

Dental Technician: Include in name

Dentist: Include in name

Patient: Include in name

Notes: Include in name

Patient ID: 120116_093910

Include time stamp in name

Standard

OK Abbruch

Figure 31: "New project"

Information such as laboratory, dental technician, patient, etc. is entered in these entry fields. The information entered here is added by setting the checkmark in the "New project" of the patient ID (Project name). If you do not enter anything in the entry fields, the time and date data will be used as standard.

After confirming with "OK", a new field will appear, allowing you to name the following measurement.

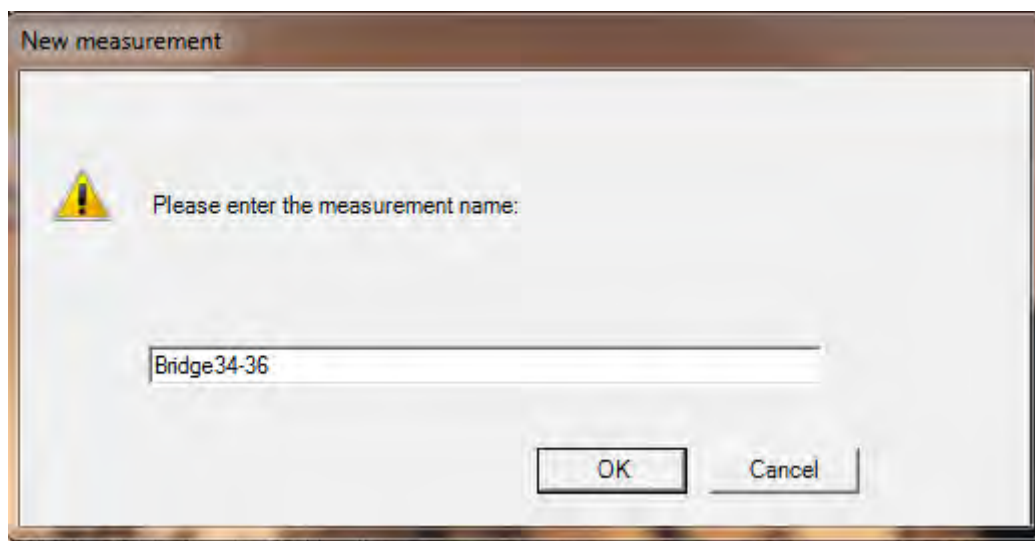




Figure 32: "New measurement"

After providing a name for the measurement, confirm this with "OK".

The scanner then begins with a preview image. After a few seconds, a 2D scan will appear on the viewer, and you can select your individual scan type.

	During the prescan or the 3D scan, the  symbol is displayed, and you can use this to cancel the corresponding scanning process.
---	--

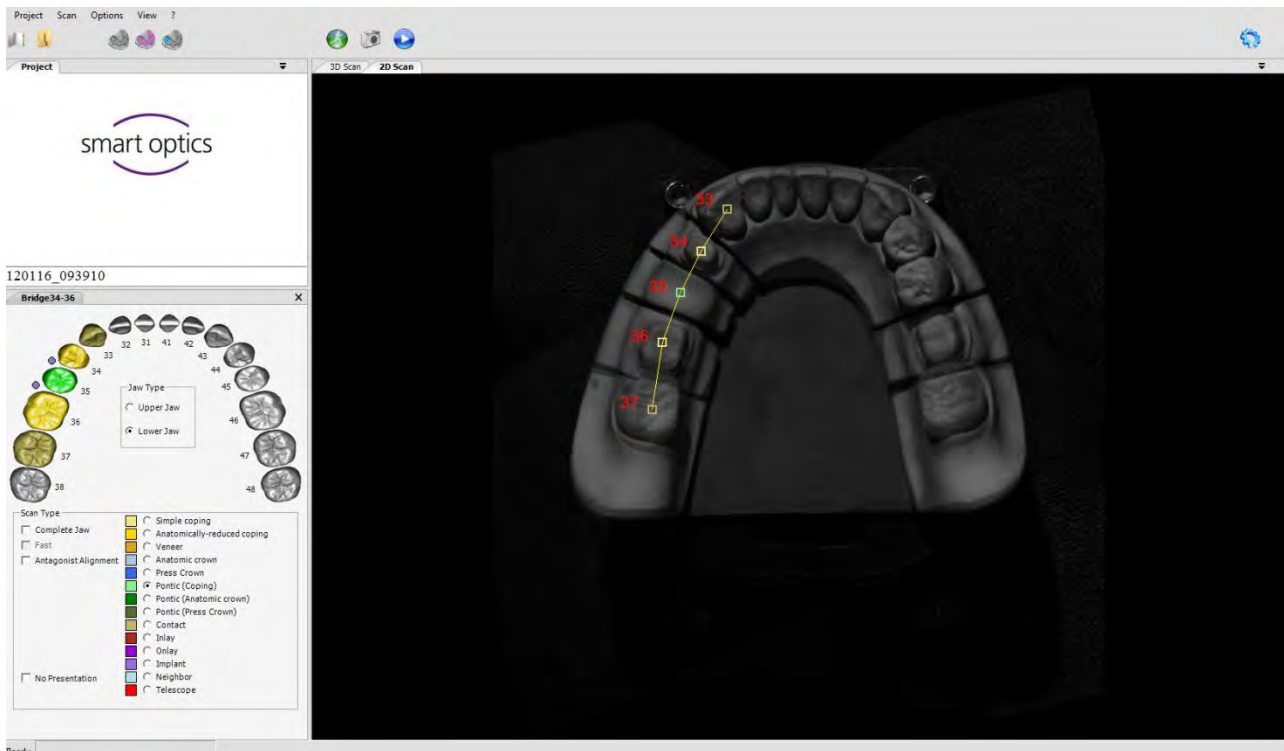


Figure 33: 2D scan and entry screen

A tooth diagram of the jaw to be scanned is located at the upper left. Activate the "Jaw type" checkbox to select whether an upper or lower jaw is involved.

Indicate under "Scan type" whether only one jaw model (standard) is to be processed or whether you have more information to be scanned, such as a bite registration, a wax-up or a scan body (abutment).

Below the tooth diagram, you will see several boxes that are coloured differently to make them easier to distinguish. These checkboxes are used to activate scan strategies for the restoration types, such as "Simple coping", "Anatomically reduced coping", „"Veneer", etc. For this, please first select your required restoration types and then click on the corresponding tooth in your tooth diagram.

You can learn about further procedures for operating the device based on the examples below.

Case study:

In this example we shall scan a **bridge of 34-36**.

Select "**Anatomically reduced coping**" under scan type and click on teeth 34, 36 and 35 as "**Reduced pontic**" in the tooth diagram. The teeth are now coloured with the relevant colour of the scan type in the tooth diagram.

Coloured squares appear in the 3D viewer at the same time (see Fig. 33).

These squares indicate the scan position of the tooth. The corresponding colouring coincides with the selected scan strategy/scan type under the tooth diagram on the left. Keep the left mouse button pressed to move these squares over the centre of the corresponding tooth stump. The scan software now detects the position of the tooth and will apply the corresponding scan strategy for these coordinates.

If you wish to delete a scan position, click on the corresponding checkpoint in the tooth diagram.

You can change the scan type of a position in the same way. Select the new scan type and click on the position that is to be changed in the tooth diagram. The scan type defined previously is then replaced by the newly selected one. The colouring of the tooth highlighted in the diagram changes correspondingly.

After aligning all scan positions, start the 3D scan with this icon.



The scanner now moves the model down automatically, based on the predefined scan strategy, and compiles a number of individual measurements from various perspectives.

In the case of adjacent teeth, it is necessary to additionally measure the individual teeth independently, as otherwise the interdental area cannot be recorded with sufficient accuracy.

The entire scene is therefore always measured in the first step, and in the second step you are prompted to present one tooth each freely.

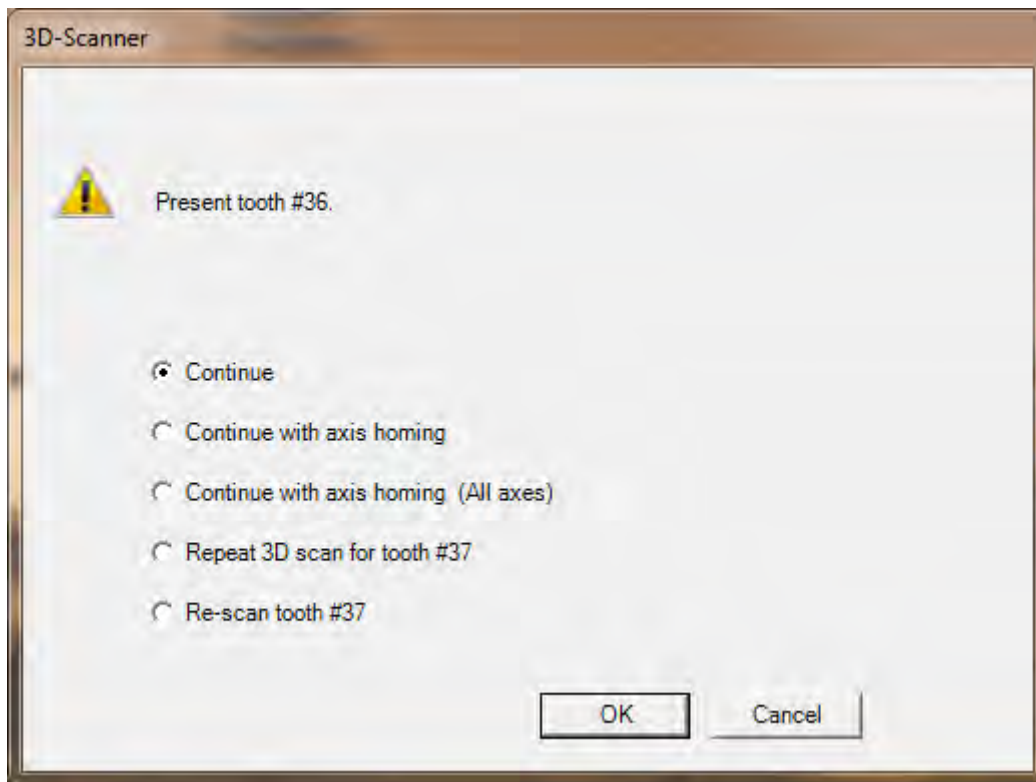


Figure 34: Prompt notification

Remove all teeth from the saw-cut model, apart from the one required.

Close the front flap and click "**Continue**" or "**Continue with axis homing**".

If "**Continue with axis homing**" is selected, a new reference travel of the axes, object holder or all axes is first executed before the measurement is continued, after removing the stumps. The object holder is moved into the predefined initial position for this.

This selection option is due to the fact that the motor driving the object holder could also be turned by hand when switched on with a correspondingly high action of force. The motor is not damaged by this, but the defined motor position is lost.

It may be the case in individual situations that greater force has to be applied for very tightly-fit stumps, and the object holder is consequently turned inadvertently. Here a considerable "**slipping through**" of the object holder can be heard, so that you usually recognize such a case. Select the option "**Continue**" to continue the scan.

After the individual tooth presentation is completed, you can evaluate the scan data in the 3D viewer for completeness (see Fig. 35).

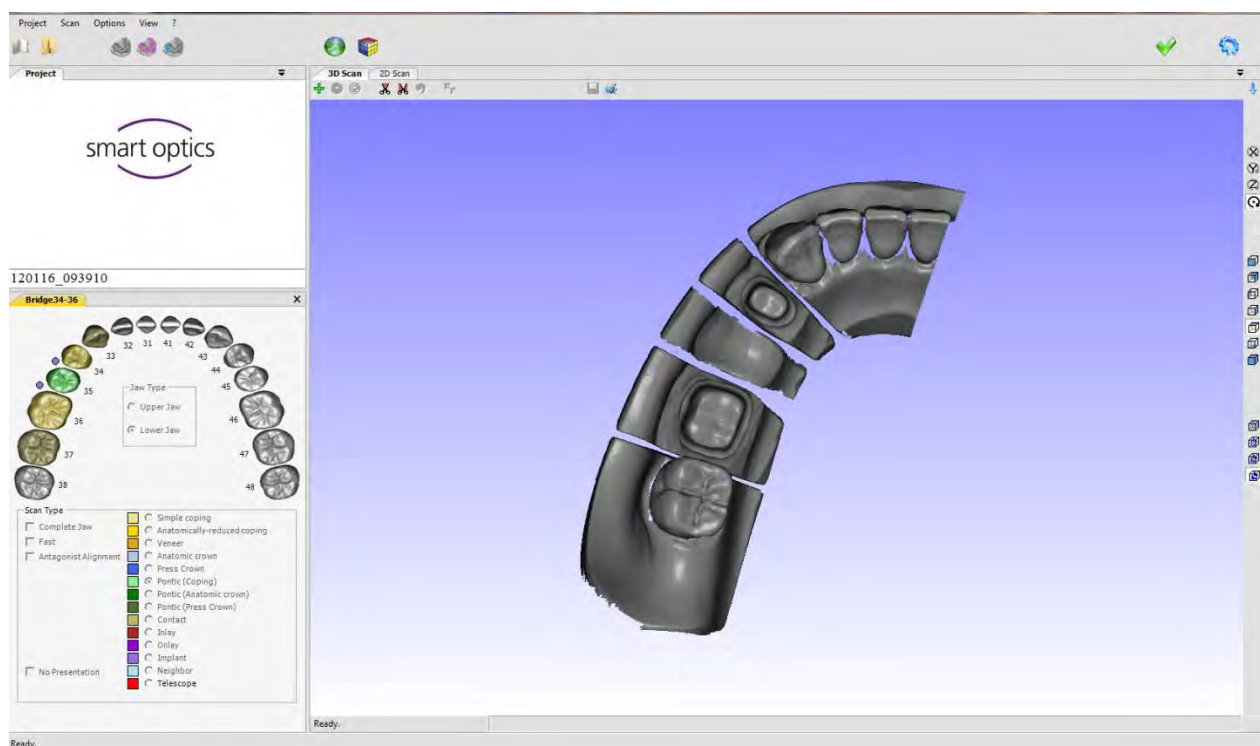



Figure 35: Data set before "matching"


	You can rotate, move and zoom the model using the following mouse functions:	
	Rotate model	Left/Right mouse button
	Move model	Both mouse buttons
	Zoom	Scroll wheel

If the model scan exhibits holes or blurred areas, you can correct these with the rescan function.


A green cross is located at the top of the 3D viewer, which can be used to start the "Rescan Mode".



If "Rescan Mode" is enabled, this green cross will be shown in the centre of the 3D viewer (see above). Now the edge of the data hole must be positioned below the cross using the familiar mouse options (left mouse button rotate, right mouse button move and scroll wheel zoom).

	If the centre of a hole is positioned directly below the cross, the scanner will move to an incorrect rescan position and possibly a recording at another point.
---	--

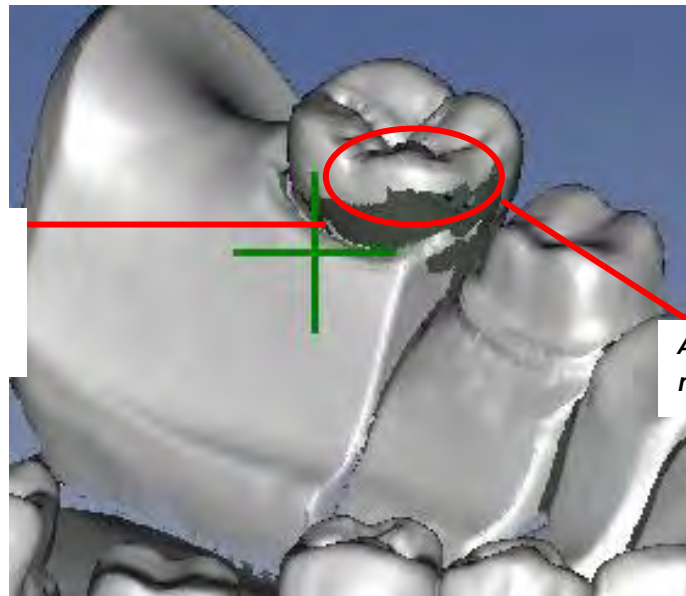Cross hairs next to the area to be scanned



Area to undergo a rescan

Figure 36: "Rescan" mode

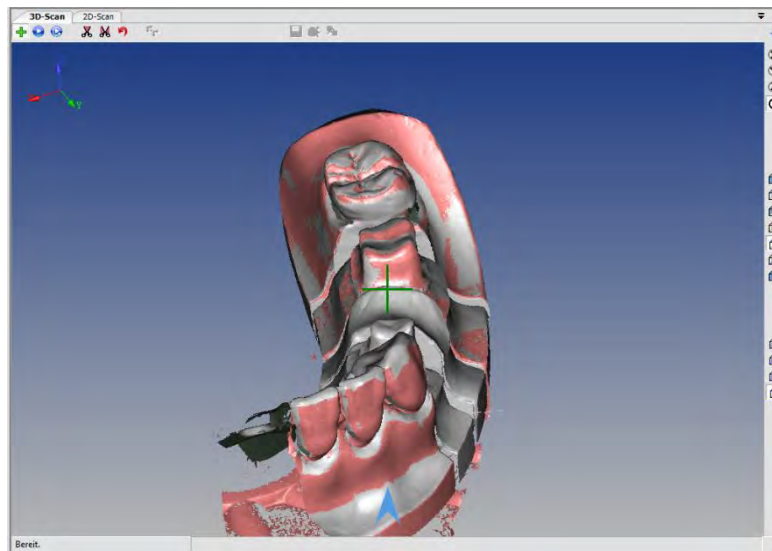

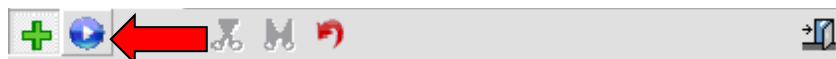


Figure 37: "Rescan"

If you click on the blue start butto  any number of rescans can now be performed until all desired areas are recorded.








For details, such as in interdental spaces, please remove all teeth that impair the sensor's view of the area to be closed, e.g. directly adjacent teeth and segments that are not required for the scan.

Once you have recorded all relevant details, click on the green cross again to complete the rescan mode.

To remove superfluous areas on the model, you can now cut the data set to size with the following tool.

-  Cut data inside the selection
-  Cut data outside the selection
-  Undo the last cutting action

Select a cut function and click with the left mouse button around the required area. In this way, you select the area on and around the model.

To complete the process, make the last selection with the right mouse button.

The model has been cut to size and can be reset with the "Undo cut action" icon.



The smaller the data set is cut to size, the smaller the data volume which has to be added together.

Save the scan data after cutting the data set to size by clicking the "Save" icon.




Figure 38: Save project

After completing the scan and cutting the data set to size, you need to perform a matching. To do this, click on the following icon:



Figure 39: "Matching" process

	<p>The recorded individual measurements are automatically added together by special algorithms, and an STL is created. This process is known as matching. You can set the quality of the STL under Item 10.4.</p>
---	--

After matching the records, the completed data set is reloaded into the 3D viewer. If more areas are missing, the rescan mode can be repeated on the STL data set.

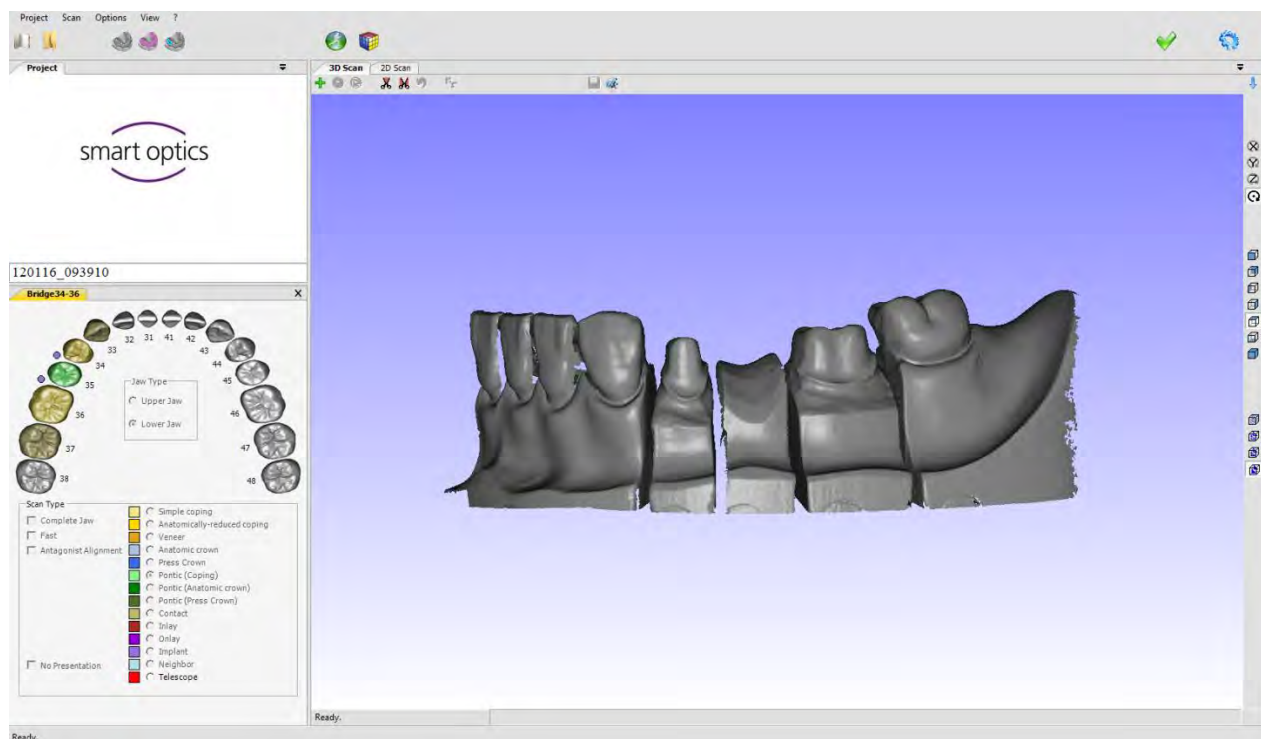


Figure 40: 3D STL data set

Once the model data set is compiled, it is possible to add a bite registration scan to the project. The bite registration function must be activated for this:



Figure 41: "bite registration" symbol

You are then prompted to position the bite registration and give the measurement a name. Once you have given the measurement a name and confirmed with "**OK**", a new information window will appear.

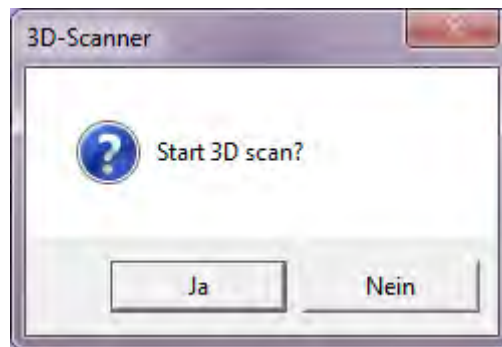


Figure 42: Information window

Confirm with "**Yes**", and the device will immediately commence with the 3D scan based on the selection created in the main scan (see Fig. 42).

After confirming the message with "**No**", the prescan is started and you can then newly define the selection for the bite registration in the prescan image (see Fig. 43)

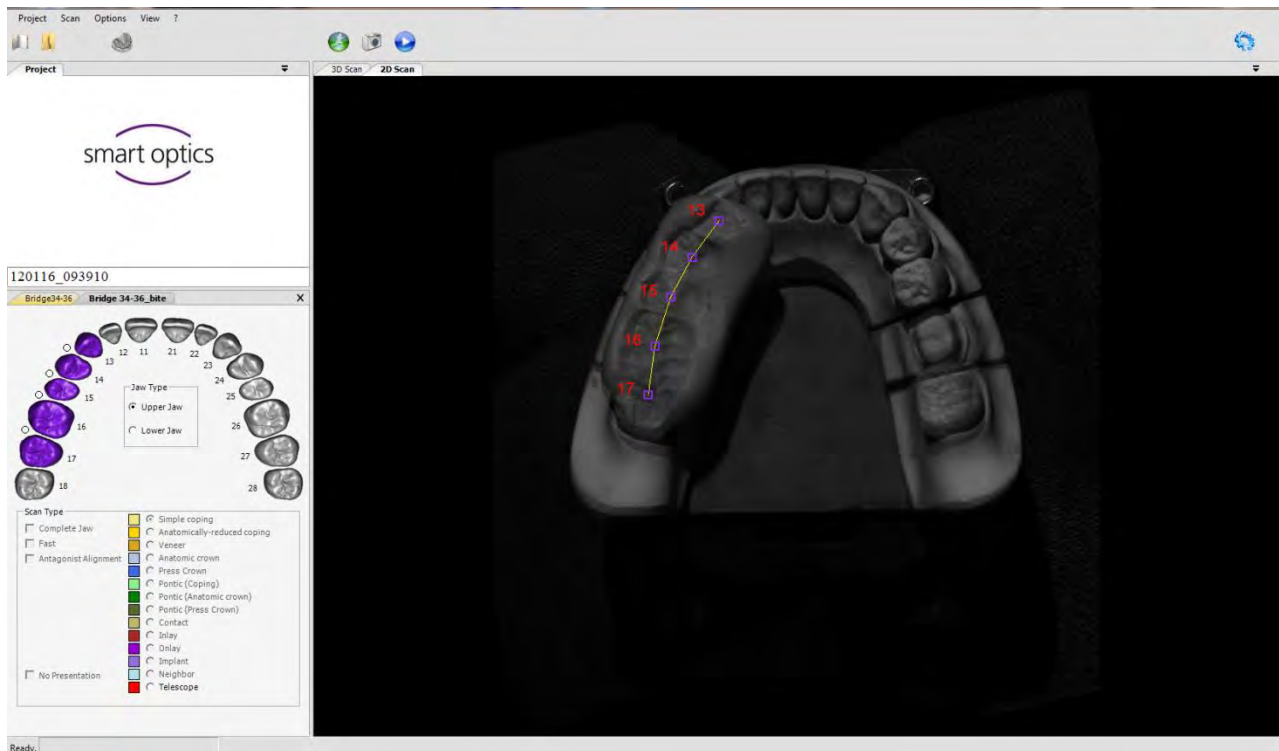


Figure 43: Squeeze bite alignment prescan



Next the 3D scan is performed as usual with the following symbol:

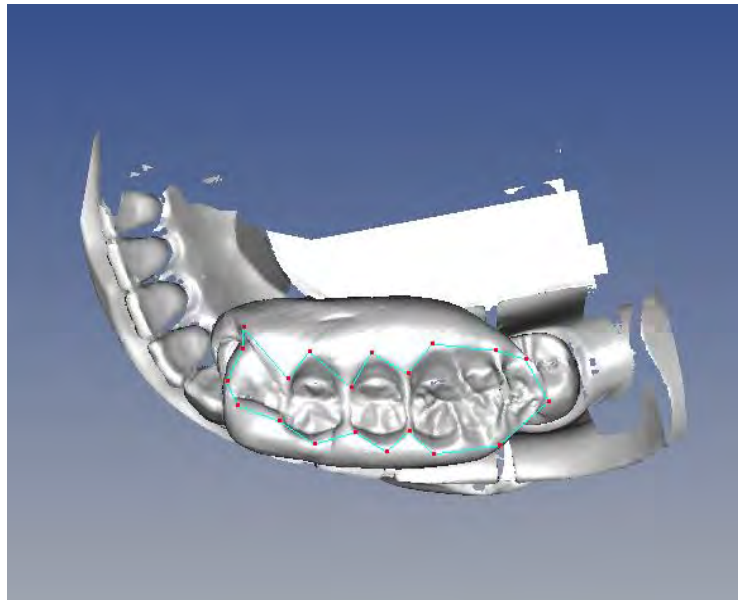


Figure 44: Cut bite registration

After completing the 3D scan, you can evaluate the result and, if necessary, perform rescans, or cut to size with the cutting tool described above. Save the change and start the matching process.



The matching button for a bite registration scan is not enabled unless the model data set has been matched previously.

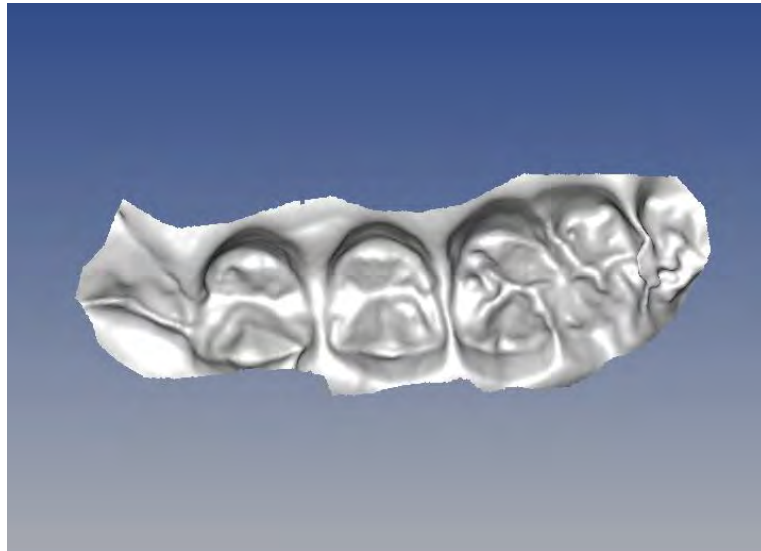


Figure 45: Cut bite registration

The bite registration is then loaded into the 3D viewer and visualized.

12.4 Case study: complete jaw scan

Clamp the model on the removable object holder. Ensure the height alignment in the measuring range. (See Fig. 30).



For certain applications, e.g. archiving in orthodontics, it is necessary to simply and quickly scan a complete jaw. Unless high detail accuracy is required for this, we recommend setting the thinning to MEDIUM or STRONG.

Create a new project as described in section 9.3 and start the 2D scan.

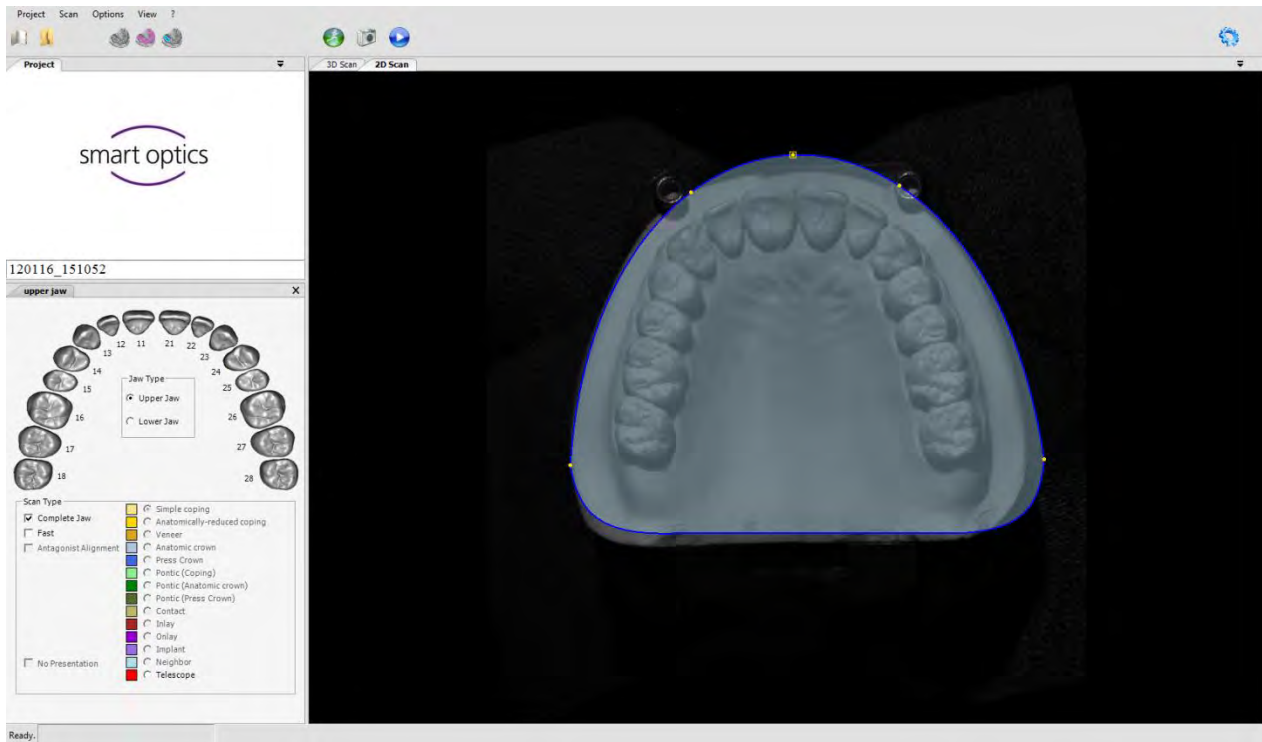


Figure 46: Complete jaw 2D scan

Select between lower and upper jaw correspondingly under "**Jaw type**".

Activate the "**Complete jaw**" checkbox under "**Scan type**" (see Figure 46).

A new control (blue line) is shown in the prescan view:

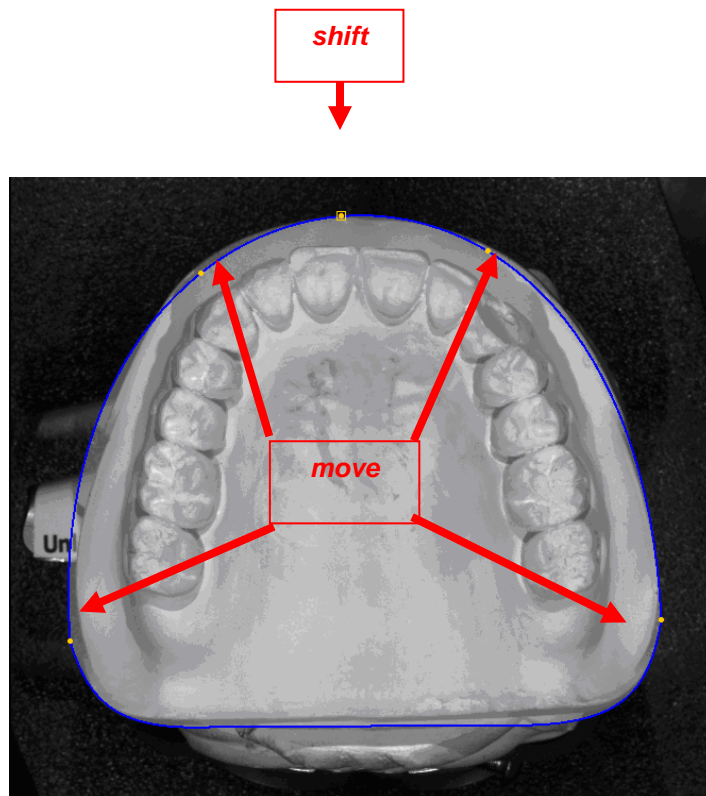


Figure 47: Complete jaw alignment prescan

The scan range is adapted to the size of the object to be scanned using the yellow checkpoints. The jaw must be inside the blue outline, the scan range.

Start the 3D scan.



Figure 48: Starting the scan process

The scanner now moves the model down automatically, based on the predefined scan positions, and compiles a number of individual measurements from various perspectives.

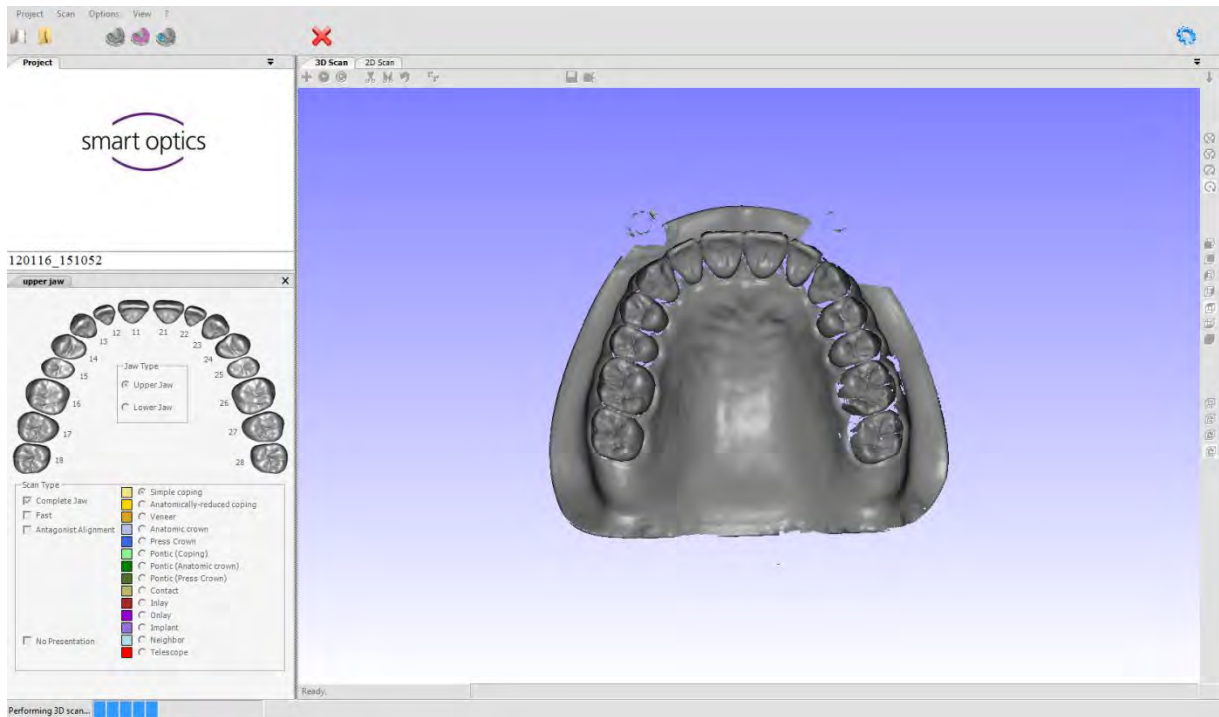


Figure 49: Complete jaw 3D scan

After the 3D scan, cut the data set if necessary. Now start the STL calculation with the Matching button.



The matching process for a complete jaw can take up to several minutes, irrespective of the number and size of the scans to be matched or the selected thinning level and computer capacity.

The jaw scan is visualized in the 3D viewer as shown below:

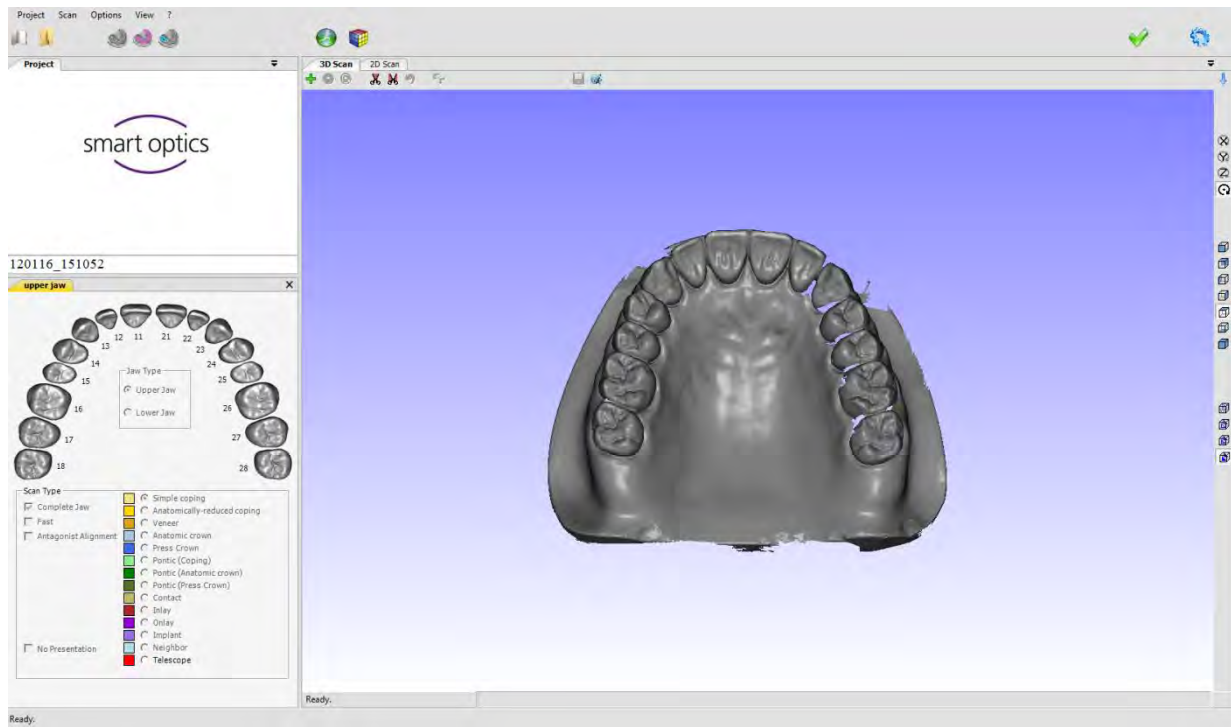


Figure 50: 3D-STL data set



You can tell from the colouring of the file tab whether a data set is already "matched", i.e. whether there is an STL. (Blue - individual records, Yellow - STL data set, Grey - 2D scan) (see Fig. 53)

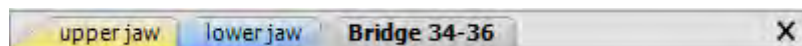


Figure 51: Tabs

12.5 Case study: occlusion scan

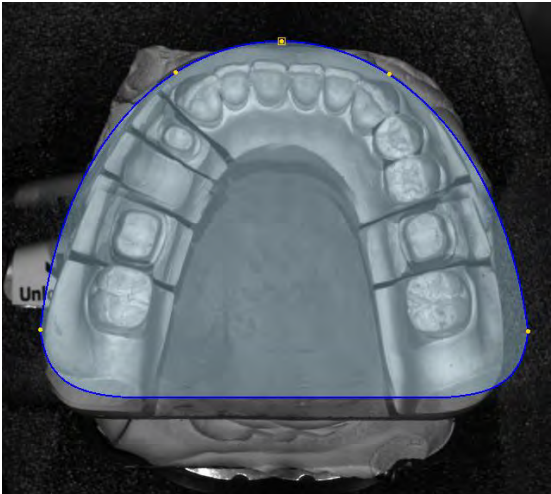
You can also scan a jaw pair in occlusion.

In our example, we shall also scan the lower jaw scan from section 9.3. (3-member bridge) and the upper jaw scan from section 9.4. (complete jaw scan) in occlusion.

To do this, we compile the scan data for the upper jaw and the lower jaw as described above. Matching is not yet carried out.

Sample images:


Prescan lower jaw



STL of the lower jaw



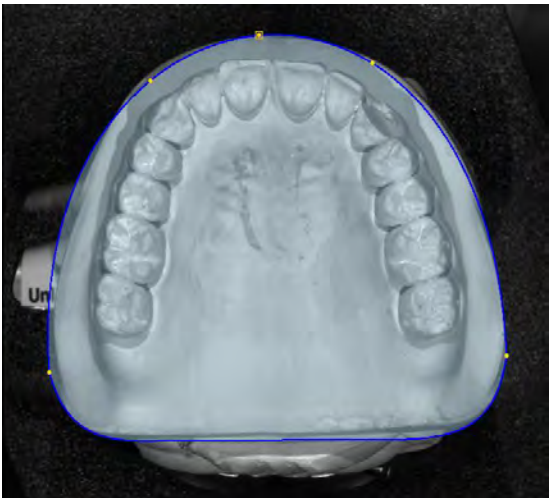Prescan upper jaw



STL of the upper jaw



As soon as you have scanned in the lower and uppers jaw, you must perform a "vestibular scan", which allows the software to calculate how the jaw halves are located relative to one another.



The scanned jaws do not have to be matched, as all scans performed are added together at the end of the process.

You can align the models with an elastic band, adhesive stick or other fixing option, so that the models sit rotationally secure to one another.

Place the models fixed to one another in the scanner and close the lid.



Use the spacer plate to ensure the correct height alignment.

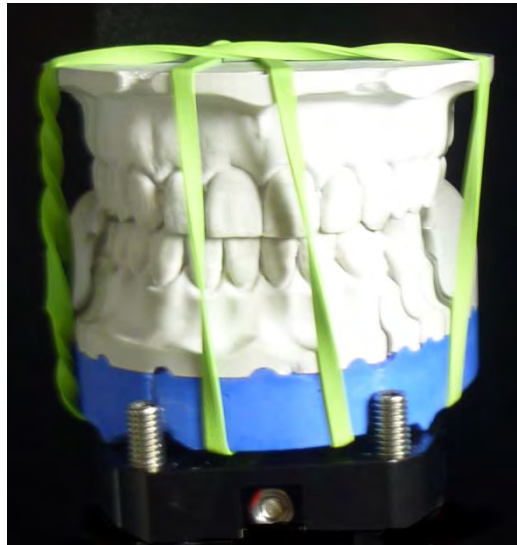


Figure 52: Fixed jaw

You must activate "**Add dental scan**" by clicking on the following icon to carry out a vestibular scan:



Figure 53: Add dental scan

Give the file a name and perform a prescan.

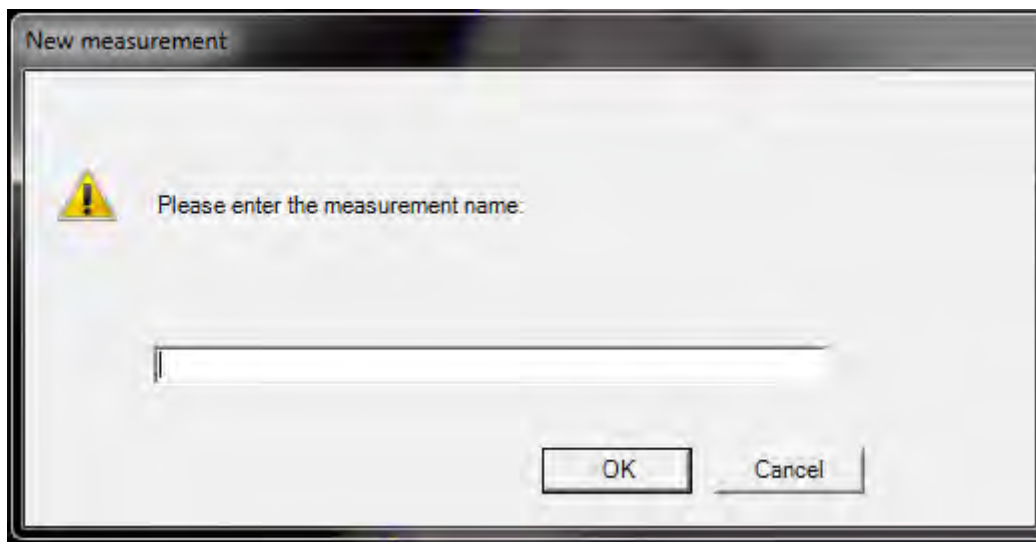


Figure 54: New measurement

Once the prescan has been performed, click on "**Antagonist Alignment**" in the left of the project entry screen, after which you can select what method is to be used to perform the vestibular scan.

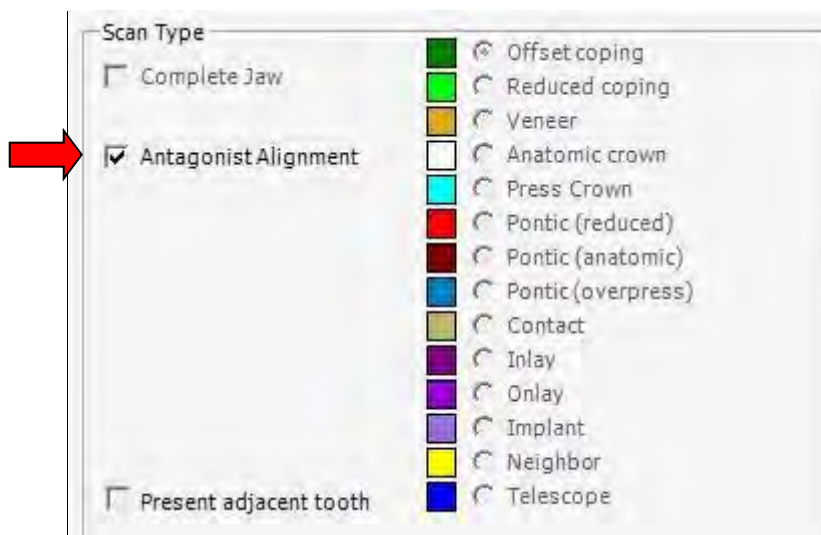


Figure 55: Antagonist Alignment

A blue outline will appear in the 2D image after clicking on "**Antagonist Alignment**".
Align this to the vestibular progression of the 2D scan (see Fig. 55).

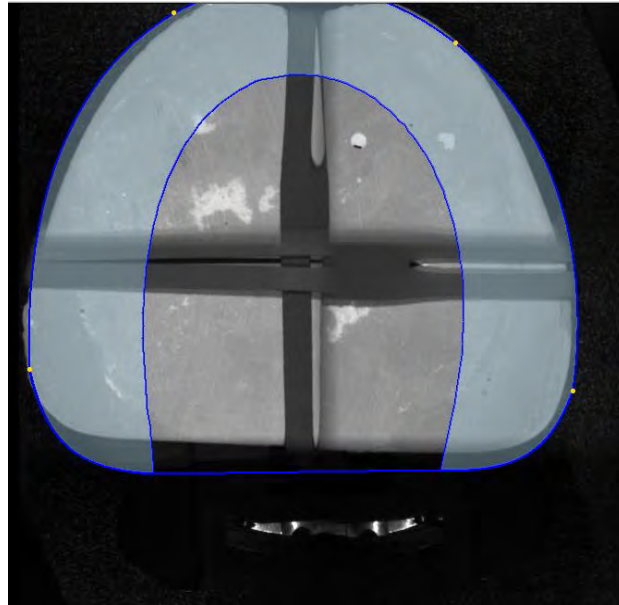


Figure 56: Prescan

Start the scanning process as usual.

The outer sides (vestibular view) are scanned by the device, as shown below:

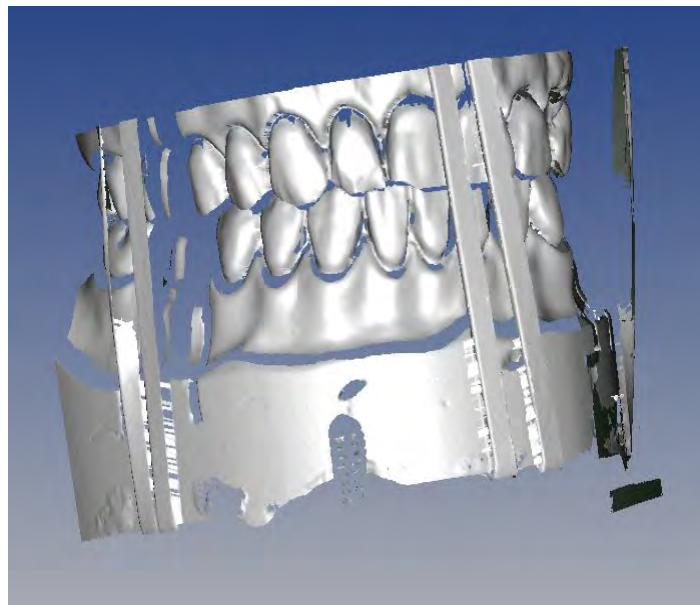


Figure 57: Vestibular scan

Click on the articulation icon:



Figure 58: "Articulation" icon

The software now indicates in succession which jaw has to be loaded.

First confirm the "Lower jaw tab" with OK. Then proceed in the same way with the "Upper jaw tab".

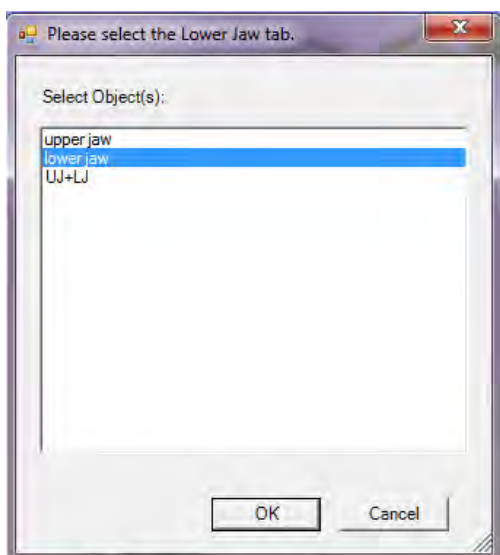


Figure 59: Select lower jaw

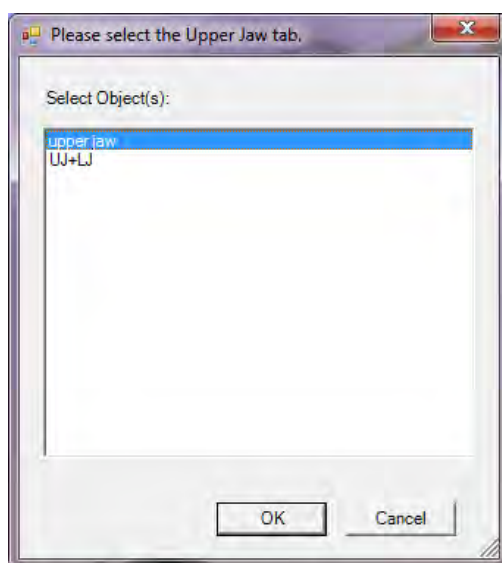


Figure 60: Select upper jaw



The alignment for the "NOARTICULATOR" must always be performed with the 3-point method.



To obtain a complete, mutually aligned STL data set from both jaw halves, the individual scans must be added together. You have to indicate reference points on the scans for this.

Select at least three reference points by clicking on them once with the left mouse button. Points which you can find again on the vestibular image are ideally suitable for this.



If you have clicked on a point by mistake, you can undo this by with the "**Delete**" key on the keyboard.

You can move the scan image with the right mouse button and rotate the view in the 3D viewer with the "Shift" key and the left mouse button.

Once you have highlighted at least three points, press the "**space bar**" on your keyboard to toggle between the lower jaw scan and vestibular scan.

Search for and select the same points as previously for the vestibular scan.

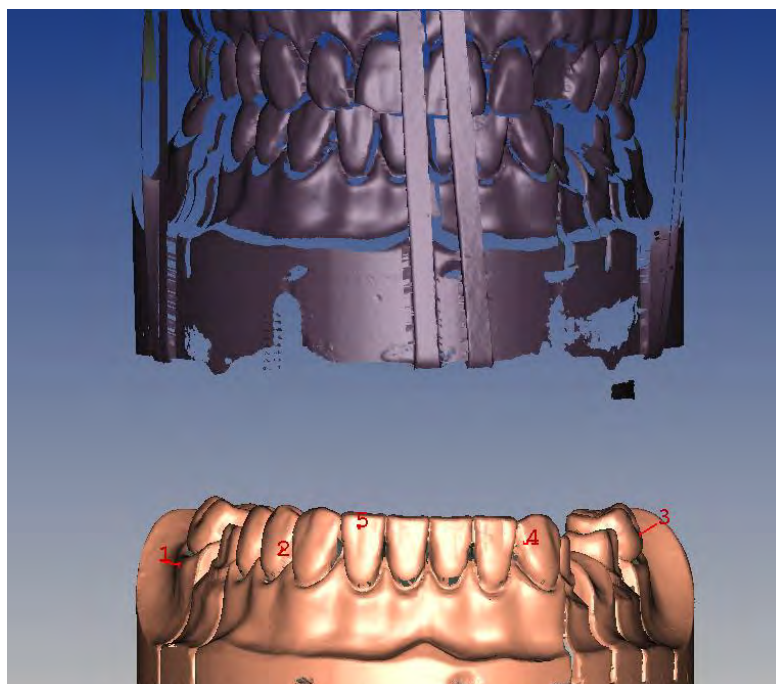


Figure 61: Lower jaw alignment

Once you have selected the reference points of the lower jaw to your satisfaction, click on the "**Next**" icon in the viewer to load the upper jaw.



Figure 62: "Next"

Follow the same procedure as with the lower jaw and select the "Next" icon again once you are satisfied with the selection.

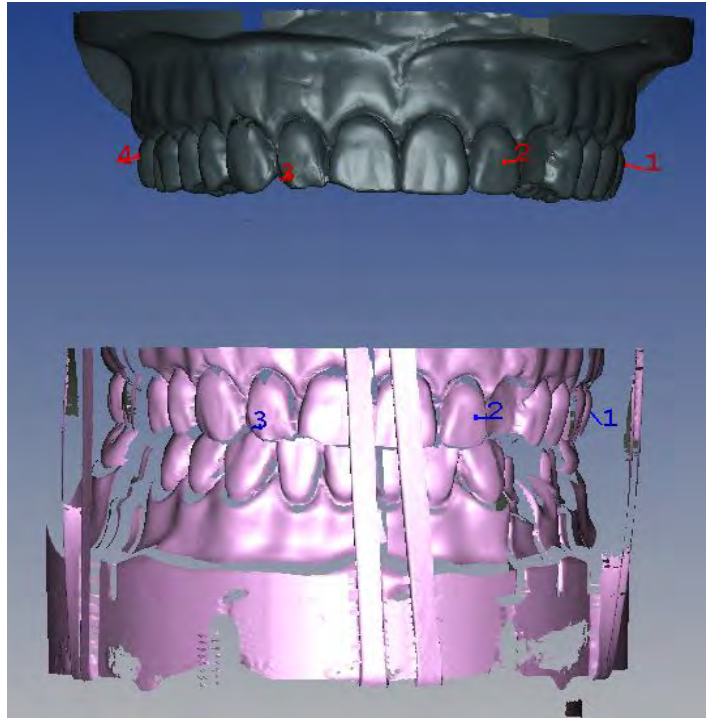


Figure 63: Upper jaw alignment

A calculation for the jaw is now carried out. After a few seconds, the result will be displayed to you in the 3D viewer, which you can then check for correctness.

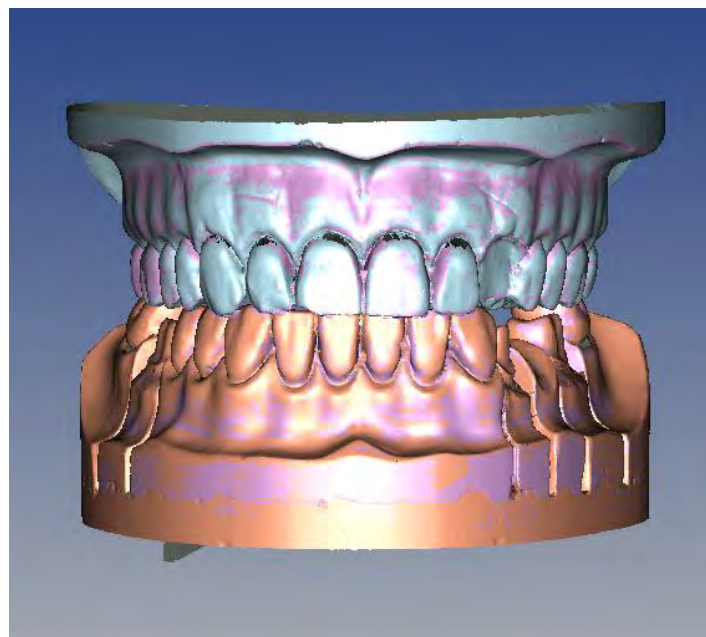


Figure 64: Aligned upper and lower jaws



You can detect an offset more easily using differing colour schemes for the upper and lower jaw.

If you are pleased with the data set, confirm with the "**Next**" button, after which the software will automatically commence the matching process.

All data sets performed are now matched. This can take several minutes.

Sample image for a finished STL data set in an upper versus lower jaw situation:

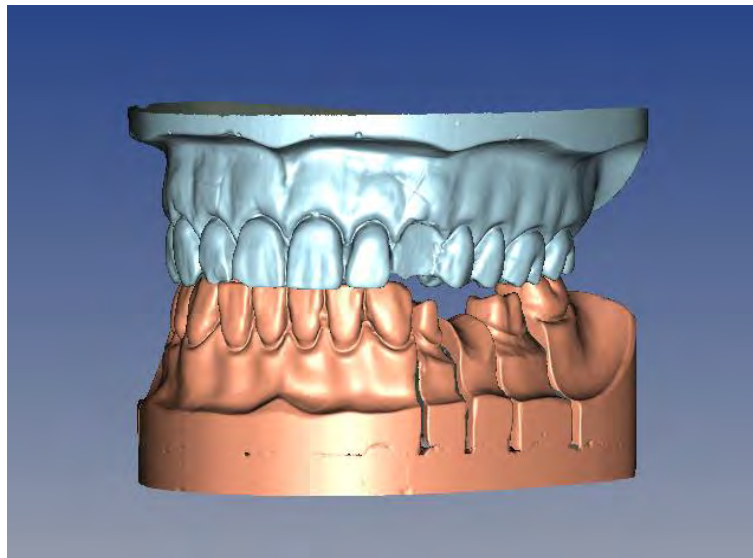


Figure 65: Articulated STL data set



The data sets can now be loaded into various CAD programs for further processing.

12. Symbols

13.1 The symbols



Activity icon



Perform 3D scan



Perform 2D scan



Move automatically to the service position



Open a new project



Open an existing project or STL file



Start the matching process



Add a new measurement



Add bite registration



Add wax-up



Continue in workflow



Cut data inside the selection



Cut data outside the selection



Undo the last cutting action or measurement



Open the dialog "Fill holes"



Save the last work step (possible if activated blue)



Open "Settings" dialog



Information on the product



End the Activity program

13.2 The 3D viewer symbols (object view)



Show camera angle



Show front view



Show rear view



Show left side



Show right side



Show upper side



Show lower side





Show data set in isometric alignment



Rotate the object around the X axis only



Rotate the object around the Y axis only

	Rotate the object around the Z axis only
	Normal rotation mode

13.3 Activity Menu options

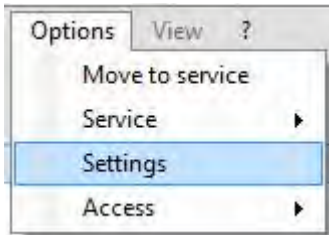


Figure 66: Activity Menu options

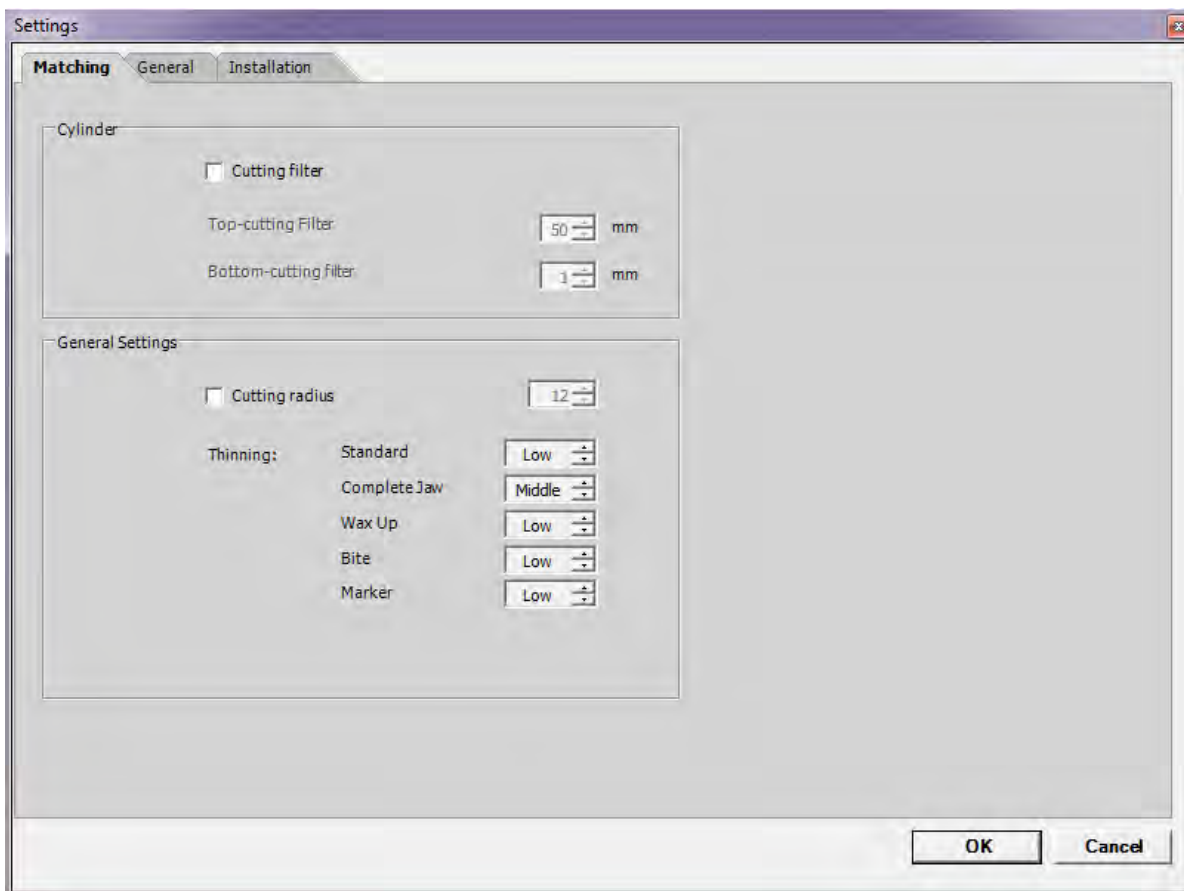


Figure 67: Settings - Matching

Cylinder:

The checkbox is used to activate a **cutting filter**, which cuts the upper and/or lower area of the STL data set. The dimensional unit is millimetres.

General settings:

Thinning out is used to determine the STL accuracy and hence the data size.

The maximum resolution of the sensor is reproduced at "OFF".

A somewhat larger spacing of the pixels to a data set is linked under "LOW".

This spacing is enlarged further at "MIDDLE" and "STRONG".

"AUTO" varies the thinning out as selected. The stronger the thinning out, the smaller the linked data set. The greatest accuracy is therefore attained when thinning out is set to "OFF" or "LOW".

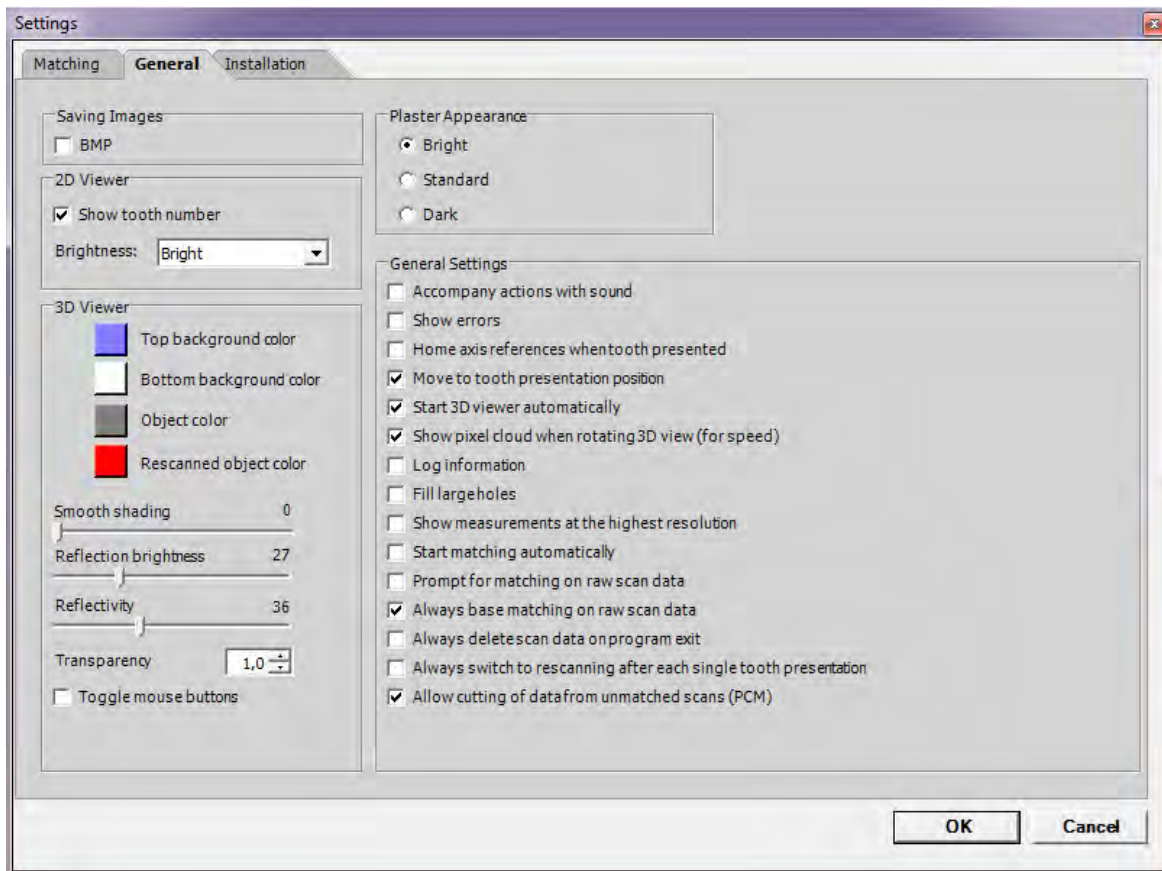


Figure 68: Settings - General

Saving images

BMP:

A **BMP** (bitmap) of the relevant recording position is created upon activation. This image, in which the light stripe pattern can also be seen, is also saved in the project directory. These images serve as an aid if any measuring errors occur.

2D viewer

Brightness:

You can set the light intensity of the camera live image here. These settings do not have any effect on the 3D measurements.

3D viewer

Colour top:

You can set the background colour of the viewer for the upper area here.

Colour bottom:

You can set the background colour of the viewer for the lower area here.

Object colour

You can determine the colour in which the 3D object is to be displayed here.

Object colour rescan

You can determine the colour in which the last manually added image is to be displayed here.

Smooth shading

Smooths the surface of the 3D object in the viewer additionally. This setting does not have any effect on the 3D measurement.

Reflection brightness

You can adapt the reflection of the light on the object in the 3D viewer here.

Reflectivity

You can determine the intensity of the reflection on the object in the 3D viewer here.

Transparency

If you activate the checkbox, you can change the button command assignment on your mouse for moving and rotating the object in the 3D viewer.

Plaster appearance

The light intensity of the sensor during the measurement is influenced by this. "Bright" for white plaster, "Standard" for beige-coloured, and "Dark" for very dark plaster.

General settings:**Accompany actions with sound**

If activated, the computer will emit a signal for every prompt via the installed speaker. However, this function is not supported by every PC system.

Show errors

If you activate "show errors", a window with the corresponding error message and description will appear if exceptions occur.

Show pixel cloud when rotating 3D view (for speed)

If you activate the checkbox, the object will be displayed in the 3D viewer as a pixel cloud when moving with the mouse. This optimizes the visualization speed during the movement.

Log information

If activated, every work step of the scanner and software as well as the joint communication is saved in a special folder.

This function is necessary for technical support.

A log file is written that logs all settings and procedures and issues a detailed error message if problems occur.

Fill large holes

If sporadic "data holes" appear after a scan, you can close these using the "Fill holes" function. The size of the hole to be closed is limited upwards by entering the area in mm². We recommend that you use this function for flat surfaces only. Incomplete areas at preparation limits or sharp, occlusal edges should not be extended by this, as the area to be replaced is only interpolated.

Show measurements at the highest resolution

If you activate the checkbox, the object will be displayed in the 3D viewer at an even higher resolution. The object calculation when adding new records as well as during movements requires much greater capacity from the graphics card and can lead to time delays.

Prompt for matching on raw scan data

You can use this function to define whether matching is to be made to the existing STL

or to the existing individual images. The query is only set for an existing STL.

Always base matching on raw scan data

To create the STL, the existing individual images are used again during rescans or supplementary scans, and an existing STL is never added for calculation.

Always delete scan data on program exit

The existing individual images are automatically deleted irrevocably when the program closes.

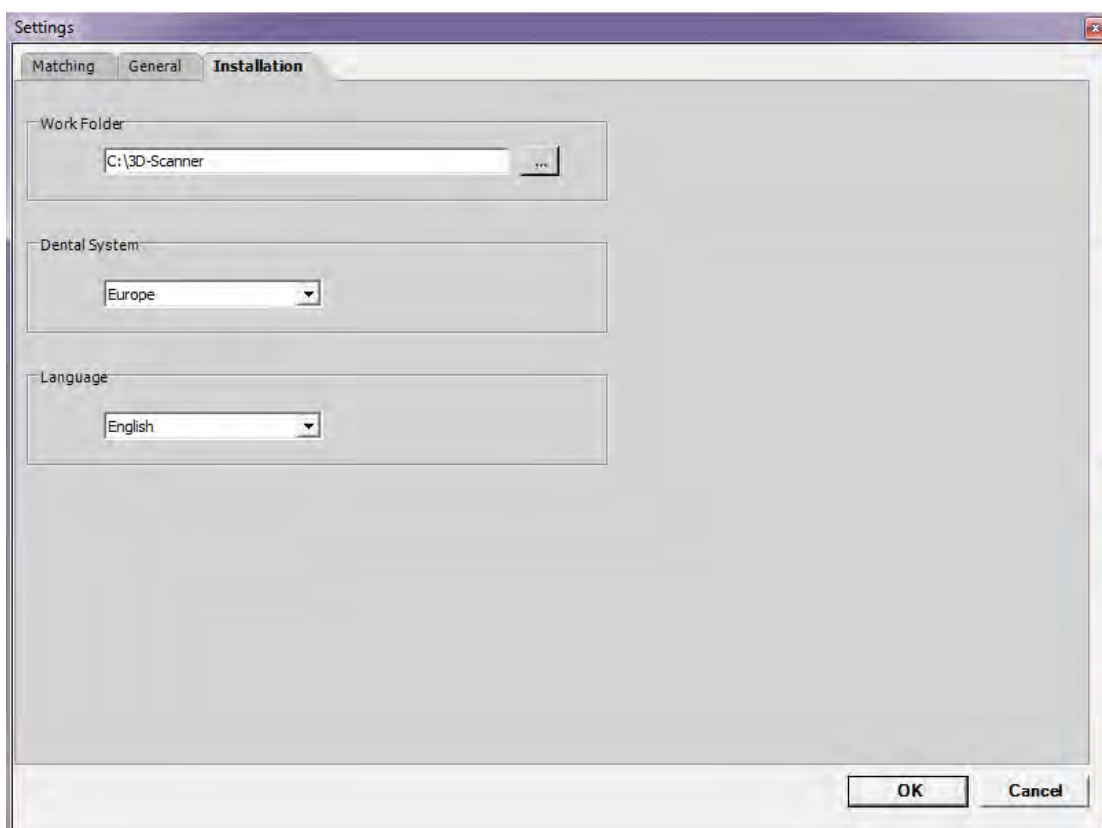


Figure 69: Settings - Installation

Work Folder:

The storage location of the scan data can be changed here. All data is created and saved in this folder.

Dental System:

You can indicate what tooth scheme you prefer in the Dental System.

Language:

The language of the user interface is set here. The following are available for selection:

German, English, Spanish, French, Italian, Portuguese, Romanian, Turkish, Chinese, Chinese Traditional, Greek and Russian.

13. Maintenance and servicing

The axes should be calibrated every four weeks using the calibration model supplied to ensure consistently good results. They should also be calibrated each time the device is transported. We also recommend performing an axis calibration at temperature fluctuations of +/- 15 °C.

The scanner should be cleaned regularly when in service. To do this, turn the scanner off and use a vacuum cleaner to carefully remove any dust and foreign objects that have accumulated in the base area of the scanner.

The optics of the 3D sensor are situated in the upper area of the interior of the scanner. Do not attempt to clean them, as inappropriate cleaning may result in damage.

Cleaning with a slightly damp microfibre cloth is sufficient in most cases. Do not use paper towels, etc., if at all possible, as they can easily scratch the sensitive plastic surfaces.

Please do not treat the case with strong cleaning agents.

14. Faults and repairs

Please remember that the scanner is a sensitive optical device. Maintenance and repair may therefore be performed only by trained technicians.

Faults that cannot be cleared by restarting the scanner or the program should be referred to customer service.

15. Environmental factors and disposal

16.1 Packaging

You can return the packaging to your dealer for disposal.

We strongly recommend that you keep the packaging so that it can be used if you need to transport the scanner or if you have to return the scanner to us in the event of a warranty claim.

16.2 Disposal

You can return the scanner to the manufacturer or dealer for disposal.

Please remember that the scanner is designed exclusively for commercial or industrial use.

It must therefore not be taken to public waste facilities for disposal.

Please contact your dealer or the manufacturer regarding disposal.

You can find out more at www.smartoptics.de by going to Company and then Environment.

WEEE registration number: DE47893210

16. Imprint

Manufacturer:



smart optics Sensortechnik GmbH

Lise-Meitner-Allee 10

44801 Bochum

Germany

Please contact your dealer if you have any questions.

We reserve the right to make changes due to technical progress.