IB-8080

Digital Breast Beauty Instrument INSTRUCTION



Dear customer,

Thanks for your choosing our product!

For the best capability and safety, please read the instruction before operation. Please reserve the instruction book for thereference.

Principle and effect

Principle of micro electricity

When micro electricity passes through body, it can form a magnetic field. When electromagnetic wave goes through breasts, it makes a kind of tingle feeling, that is because the magnetic field is dredging the block part inside body, activate gene of mammary glands and make it re-grow, get through lymphatic, accelerate lymph circulation, in this way to rebirth breast fat and get the effect of breast enlargement.

• Principle of rhythm physical suction

- 1. It employs bionic technology, lead estrogen to breast and make the breast to be action target, in this way to make breast fat accumulate and enlarge the breasts.
- 2. The strong movement can give muscle a plentiful exercise, meanwhile stimulate the perspire of estrogen, in this way to get a natural breast care.
- 3. The suction can give exercise to musculature and improve its elasticity, so that to lift breast and make it firm. The suction and releasing movement can give stimulation to chest muscle and mammary glands, finally stimulate the growth of breast.
- 4. The function can get through lymphatic, accelerate lymph circulation, which can get a good effect to eliminate tumefaction, defend breast disease, balance incretion.

Principle of vibration

Micro wave vibration can give a continuous activation to base of mammary glands, make sport to undergrowth breast or saggy breast and make it re-grow.

Principle and Notice

• Effect

Vibration massage: Accelerate the circulation of blood and lymph, give exercise to muscle, relax muscle, easy tire.

Photon: Make skin tender and smooth, promote the essence absorbed by body

Breast care: Enlarge breast, lift breast, modify nipple, modify chest, recover breast elasticity

Face care: Thin face, suck blackhead, reduce dark eye-socket, remove eye pouch, modify double chin, lift saggy canthus.


Back care: Vacuum metabolism, skin-scratch treatment, acupressure therapy, lymph de-toxin, frap skin, strengthen skin flexibility

Obesity cure: lose weight, draw in abdomen, thin waist

Notice:

the straw can not suck water. After clean the connector and the accessories please make sure there is no water inside before using.


Introduction of panel and accessory

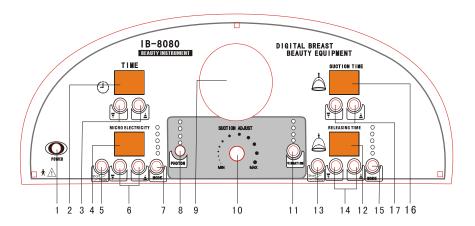

- 1. Power
- 2. Time display
- 3. Time adjust (0-60minutes)
- 4. Micro electricity intensity display
- 5. Micro electricity start/pause
- 6. Micro electricity intensity adjust (intensity 1-9, from weak to strong)
- 7. Micro electricity mode select
 - Sine wave mode
 - Cell activation mode
 - Pulse wave mode
 - Movement mode
 - Surface massage mode
- 8. Photon mode select
 - Random light
 - Sparkle interval 0.5 seconds
 - Sparkle interval0.2 seconds
 - Constant light

11. Vibration mode select

- O—Pulse vibration
- Vibration 2 seconds, stop 0.5 seconds, interval vibration
- Vibration 1 second, stop 1 second, interval vibration

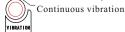

- 12. Releasing time display
- 13. Patterfunction start/pause
- 14. Releasing time adjust
- 15. Patter mode select
 - O—Automatic recycle mode2 (for breast modify)
 O—Automatic recycle mode 1 (for breast lift)
 - Suction 3 seconds, releasing 1.5 seconds(for breast enlargement)
 Suction 1 second, releasing 1 second(for relaxing massage)
 - Continuous suction(for cupping)
- 16. Suction time display
- 17. Suction time adjust


Suction adjust instruction

- 9. Suction display 2. Salast continuous
 - 2. Select continuous suction mode
- 10. Suction adjust 3. Connect pipe to machine, block the pipe as the picture, adjust the proper suction intensity



Introduction of panel and accessory



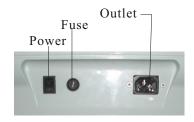

Vibration breast cup(6 pieces)






Breast cup wire Single connection pipe Double connection pipe





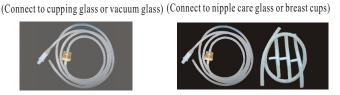


Cupping glass Nipple care glass



Connection instruction

Double connection pipe connect to breast cups

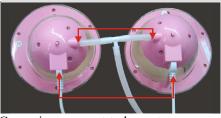



Cup wire connect to breast cup

Cup wire connect to mainframe

Double connection pipe connect to Cupping glass



Single connection pipe Pipe connect to Cupping glass Pipe connect to mainframe

Double connection pipe connect to nipple care glass




Double connection pipe Pipe connect to nipple care glass Pipe connect to mainframe

Single connection pipe connect to vacuum glass



Single connection pipe Pipe connect to vacuum glass Pipe connect mainframe

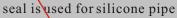






make sure that it is in the position to get full power

I. Breast care

Step 1: clean skin

Step 2: massage direction

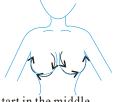


Start in the middle

from inside → down → outside

→up**→**inside

massage in a circle way

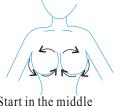

Start in the middle

from inside → up → outside

→down → inside

massage in a circle way



From outside to inside, lift the

breasts up

Step 3: Operation guide

- 1. Spread essence oil or breast cream on breast.
- 2. Paste on nipple patch
- 3. Instrument adjust
 - A. Instrument adjust for breast lifting
- (1) Select automatic recycle mode 1

-Automatic recycle mode2 (for breast modify)

-Automatic recycle mode 1(for breast lift)

Suction 3 seconds, releasing 1.5 seconds(for breast enlargement)

Suction 1 second, releasing 1 second(for relaxing massage)

Continuous suction(for cupping)

(refer to panel introduction No.15)

(2) Select sine wave mode(could be chosen)

Sine wave mode

Cell activation mode

Pulse wave mode

Movement mode

Surface massage mode

(Refer to panel introduction No.7)

(3) Select frequent vibration

- Pulse vibration
 Vibration 2 seconds, stop 0.5 seconds, interval vibration
 Vibration 1 second, stop 1 second, interval vibration
 Continuous vibration
 (Refer to panel introduction No. 11)
- (4) Select random light
 - Random light
 Sparkle interval 0.5 seconds
 - Sparkle interval0.2 seconds
 Constant light

(5) Suction intensity adjust to -20



(Refer to panel introduction No.8)

- B. Instrument adjust(for breast enlargement)
- (1) Select suction 3 seconds, releasing 1.5 seconds,
 - Automatic recycle mode 2 (for breast modify)

 Automatic recycle mode 1 (for breast lift)

 Suction 3 seconds, releasing 1.5 seconds(for breast enlargement)

 Suction 1 second, releasing 1 second(for relaxing massage)

 Continuous suction(for cupping)

(Refer to panel introduction No.15)

(2) Select cell activation mode

Sine wave mode
Cell activation mode
Pulse wave mode
Movement mode
Surface massage mode

(Refer to panel introduction No.7)

- (3) Select continuous vibration
 - Pulse vibration

 Vibration 2 seconds ston 0.5 seconds interval

Vibration 2 seconds, stop 0.5 seconds, interval vibration
Vibration 1 second stop 1 second, interval vibration

Continuous vibration

(Refer to panel introduction No.11)

(4) Select constant light

- Random light
- Sparkle interval 0.5 seconds
 - Sparkle interval0.2 seconds Constant light

(Refer to panel introduction No.8)

C. Instrument adjust (for breast modify)

(1) Select automatic recycle mode2

- -Automatic recycle mode2 (for breast modify)
- -Automatic recycle mode 1(for breast lift)
- Suction 3 seconds, releasing 1.5 seconds(for breast enlargement)
- Suction 1 second, releasing 1 second(for relaxing massage)

(Refer to panel introduction No.15) Continuous suction(for cupping)

(2) Select movement mode

- Sine wave mode
- Cell activation mode
- Pulse wave mode
- Movement mode

Surface massage mode

(Refer to panel introduction No.7)

(3) Select sparkle interval 0.2 seconds

- -Random light
- Sparkle interval 0.5 seconds

(Refer to panel introduction No.8) Sparkle interval0.2 seconds

Constant light PHOTON

(4) Select vibration 2 seconds, stop 0.5 seconds

- Pulse vibration
- Vibration 2 seconds, stop 0.5 seconds, interval vibration
- Vibration 1 second, stop 1 second, interval vibration

(Refer to panel introduction No.11) Continuous vibration

(5) Suction intensity adjust to -10



D. Instrument adjust (for breast massage)

(1) Select suction 1 seconds, releasing 1 seconds

—Automatic wave 2 (for breast modify)

-Automatic wave 1(forbreast lift)

Suction 3 seconds, releasing 1.5 seconds

Suction 1 seconds, releasing 1 seconds

Continuous suction(for cupping)

(Refer to panel introduction No.15)


(2) Select vibration 1 seconds, stop 1 seconds

Pulse vibration


Vibration 2 seconds, stop 0.5 seconds, interval vibration

Vibration 1 seconds, stop 1 seconds, interval vibration

Continuous vibration

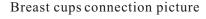

(5) Suction intensity adjust to -10

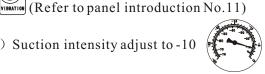

The treatment time should be adjusted to 15 minutes, adjust the suction from weak to strong, let client feel and select the proper suction intensity.

E. Instrument adjust (for nipple care)

(1) Adjust there is no light on the mode, then the parameter can be adjust

-Automatic recycle mode2 (for breast modify)

-Automatic recycle mode 1(for breast lift)

Suction 3 seconds, releasing 1.5 seconds(for breast enlargement)

Suction 1 seconds, releasing 1 seconds(for relaxing massage)

Continuous suction(for cupping)

(2) Adjust the suction time to 2 seconds, releasing time to 2 seconds



(refer to panel introduction No.16 and 17)



(refer to panel introduction No.12 and14)

(3) The treatment time should be in 2 minutes



(4)Adjust the suction intensity to-10



- (5) After finishing the treatment, clean body parts
- (6) Power off the machine, clean and sterilize the accessory



Nipple care picture

Notice: The above instrument adjust is for reference, the output parameter should be adjusted properly according to specific practice.



II. Face care

Turn on, the machine, there is only suction, no releasing, then it can make treatment of vacuum to suck comedones and flow lymph on the face.

Operation procedure

1.Vacuum

- (1). Connect the sterilized glass tube to the pipe
- (2). Start the machine and feel the suction intensity on the back of hand, adjust the intensity according to specific demand.
- (3). There are three ways to make vacuum: discontinuous vacuum, continuous vacuum, strong vacuum. When make treatment to the oily skin, hard skin, need use much stronger vacuum. When make

Treatment to xeransis skin, sensitive skin, senile skin, need much lighter

Vacuum.

- a. Discontinuous vacuum. During operation, the middle finger should cover and release on the small vent hole of glass tube in a frequent and rhythm way. The vacuum should be not too strong. This kind of operation is better to treat tender skin, loose skin and thin skin.
- b.Continuous vacuum. The middle finger cover the small vent hole of glass tube, and move together with the tube for a while then release the finger. The vacuum need much stronger. This kind of operation is better to treat oily skin, hard and thick skin.
- c.Strong vacuum. The middle finger cover the vent hole of glass tube and don't release, put the tube mouth on the skin, suck and release on the skin frequently. The vacuum intensity need strong. This kind of operation is better to treat oily skin, such parts as nose tip, nose side and the parts with comedones.

 Page 11

2. Lymph flow

Step 1: forehead

Guide: From central forehead to

the side, from up to down.

Step 2: temple

Guide: From hair side to the

bottom of ear.

Step 3: cheek

Guide: from nose wing to the

bottom of ear.

Step4: lip

Guide: from philtrum to the

bottom of ear.

Step 5: chin

Guide: fron center of chin to the

bottom of ear.

Attention: During lymph flow, should adjust the suction to weak and avoid to make the face red and

Turgescence.



Step 1



Step 2



Step 3



Step 4



Step 5

III. Body care


2. The continuous suction can do cupping action.

(1)Cupping glasses






(2) Select continuous suction mode (refer to panel introduction 15)

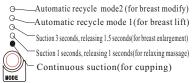


(3) Adjust time to 30 minutes (refer to panel introduction 2)





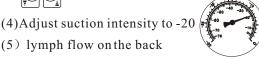

- A. Spread essence oil on the body.
- B. Put the button Start/pause to start the machine.
- C.Guide:
- C.1.from neck to ear(guide1)
- C.2.from back to oxter(guide2)
- C.3.from small back to side of abdomen(guide3)
- C.4. from arm to oxter(guide4)
- C.5. from hip to the crease(guide5)
- C.6. from backside of thigh to front side(guide6)
- C.7.from knot to front side(guide7)
- C.8.from backside of leg to knot(guide8)

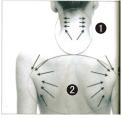



Solution for simple failure

| Failure | Solution | |
|------------------------------------|--|--|
| Dawar an main machina | Check power connection | |
| Power on, main machine no response | Check whether power source is on | |
| No output | Check the air pipe connection | |
| During operation, instrument | Using for toolong time, need to cooldown in anature way | |
| is overheat | Check whether voltage is right (need professional technician to check) | |
| No power | Check whether fuse is burnt, if so need to replace a new one | |

Technical parameter

| Electricity voltage:AC110V□ | 220V□ | 230V□ | 240V□ |
|-----------------------------|-------|-------|-------|
| Working frequency:50HZ□ | 60HZ□ | | |
| Downer 60W | | | |

Power:68W

Packing size:64×54×23CM(Table style:IB-8080)

Gross weight: 6.7 KG(Tabel style: IB-8080)

Packing size:106×55×46.5CM(With stand:IB-8080A)

Gross weight:25KG(Withstand:IB-8080A)

Packing list

| Name | Unit | Qty | Remark | |
|--------------------------|------|-----|-----------------------|--|
| Mainframe | set | 1 | | |
| 1# vibration breast cup | Pc | 2 | | |
| 2#vibration breast cup | Pc | 2 | | |
| 3#vibration breast cup | Pc | 2 | | |
| Breast cup wire | Pc | 2 | | |
| Double connection pipe | Set | 1 | | |
| Single connection pipe | Pc | 4 | | |
| Cupping glass | Pc | 2 | Big one and small one | |
| Cupping glass handle | Pc | 1 | | |
| Small nipple care glass | Pc | 2 | | |
| Middle nipple care glass | Pc | 2 | | |
| vacuum glass | Pc | 3 | | |
| vacuum glass handle | Pc | 1 | | |
| Head of tender | Pc | 4 | | |
| Pulley heads | Pc | 2 | | |
| Airproof loop | Pc | 2 | Big one and small one | |
| Power wire | Pc | 1 | | |
| Filter acrylic spanner | Pc | 1 | Stand style | |
| Fuse | Pc | 2 | | |
| Instruction book | Pc | 1 | | |

Notice

Notice: after the operation please clean the filter and head of tender by using water. If not the pump of the machine will be broken (show as picture)







Filter acrylic spanner

1. Method of cleaning the head of tender





Please hold the rubber hose Please take the head of Clean the hole with water tender away

2. Method of cleaning the filter



Lock the filter by spanner



Hold the filter metal part by hand, spin offthe plastic part by spanner



Please clean the bottle with water

Safety notice

- 1. The machine is professional instrument for beauty salon, in order to avoid any unnecessary hurt, need professional beautician to make operation.
- 2. Do not disassemble the machine or operate the machine that not explained in the instruction book. All the maintenance services shall be carried out by the professional maintenance technician approved by our company.
- 3. Do not install or operate the machine with wet hands or near the water sources so as to avoid the liquid spraying/swashing over the machine.
- 4. Pull out the plug and connect with the maintenance service center in the following cases:
 - A. liquid spray/swash over the machine or the machine get caught inrain or water.
 - B. the machine has peculiar smell, fume/smog or abnormal noise
 - C. power lines damaged or scratched.
 - D. instrument felldown or the shell damaged.
- 5. Do not put any item on the power wires, or install the instrument in the place where its power wires will be trampled.
- 6. Pull out the power wires at the end of use for the sake of safety.
- 7. Do not insert any type of item into the instrument opening; otherwise there will be a fire or electric shock. In case of something fell into the opening but you could not take it out, please pull out the power wires and connect with the dealer or our company.
- 8 .Do not put the instrument on the unstable trolley, shelf or vehicle, if fall down, the instrument may be damaged.
- 9. To guarantee the life circle of this machine, the interval of the continuous switch on and off should be over one minute.
- 10. The power shall comply with the specification specified in the label, otherwise the machine will not run orrun abnormally, even be burnt out.
- 11. The company has right to revise or amend this manual

The company reserves all the rights of interpreting the information mentioned above.

Attention:

If the supply cord is damaged, it must be replaced by manufacturer, its service agent or similarly qualified persons in order to avoid a hazard.

This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.

Children should be supervised to ensure that they do not play with the appliance.