

FOLIGAIN[®] LASER COMB

POWERED BY:
HAIRMAX[®]

FOLIGAIN.L7x[®] Advanced Laser Comb
FOLIGAIN.L12x[®] Professional Laser Comb

USER MANUAL



TABLE OF CONTENTS

Introduction	1
Classifications for Men and Women	2
Who is the best candidate? (Fitzpatrick Skin Types)	3
Warnings and Cautions	4
Risks/Testing for Sensitivity	5
How to Monitor Device Activity	6
Device Failure Mode	6
Expected Life of the FOLIGAIN® Laser Comb	6
Summary of Clinical Study in Men	7
Summary of Clinical Study in Women	8

FOLIGAIN.L7x® Advanced Laser Comb Powered by HairMax

FOLIGAIN.L7x® Advanced Laser Comb Parts	9
Getting Started	10
Replacing the Batteries	10
Proper Use	11
Treatment Instructions	12
What's included – FOLIGAIN.L7x® Advanced Laser Comb	13
Troubleshooting	14
Safety Label	15

FOLIGAIN.L12x® Professional Laser Comb Powered by HairMax

FOLIGAIN.L12x® Professional Laser Comb Parts	16
Getting Started	17
Charging Options	17
Proper Use	18
Treatment Instructions	19
What's included – FOLIGAIN.L12x® Professional Laser Comb	21
Troubleshooting	22
Safety Label	23
Power Supply Specifications	23
Cleaning Your Laser Comb	24
Electrical Safety	25
Warranty	26
Warranty Card	27

Introduction of the FOLIGAIN® Laser Comb Powered by HairMax®


The FOLIGAIN® Laser Comb is a home-use laser device that has been clinically proven to treat hereditary hair loss in men and women. The Laser Comb is an over-the-counter (OTC) device that can be sold directly to the public.

Indications

The Laser Comb models are indicated to treat Androgenetic Alopecia and promote hair growth in females who have Ludwig (Savin) I-4, II-1, II-2, or frontal patterns of hair loss and in males who have Norwood Hamilton Classifications of IIa to V and who both have Fitzpatrick Skin Types I to IV.

How the Laser Comb Works

The Laser Comb works by providing laser energy to stimulate hair follicles. For optimal results, the laser must not be blocked by the hair and must have an unobstructed path to the scalp. The Laser Comb has teeth that part the hair and allow the maximum amount of laser energy to reach the scalp.

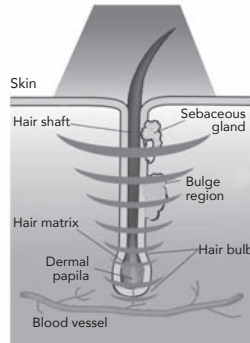

Treatment with the Laser Comb

The Laser Comb should be used 3 times per week, spaced out every other day, for example, on Monday, Wednesday, and Friday. The Laser Comb should not be used more than the recommended treatment time per day. Recommended treatment for the FOLIGAIN.L7x® Advanced Laser Comb is 15 minutes and the FOLIGAIN.L12x® Professional Laser Comb is 8 minutes. Clinical studies demonstrate that most users see benefits in as little as 16 weeks.

Laser Comb Benefits

A clinical study was conducted in men and women over a six month time period. The study showed that 93% of the men and women had an increase in the number of hairs (based on a minimum of 32 new hairs per square inch being observed at follow up).

See pages 7-8 for detailed results of the study.

It is important that you read and understand the content of this User Manual before using your FOLIGAIN® Laser Comb.


CLASSIFICATION FOR MEN AND WOMEN

WHO IS THE BEST CANDIDATE FOR THE FOLIGAIN® LASER COMB?


How do you know if The FOLIGAIN® Laser Comb is right for you?

The Laser Comb can be used by men and women with thinning hair (including colour-treated hair), or pattern baldness caused by a hereditary condition. Doctors use a system known as the Norwood Hamilton Classification (men) and the Ludwig-Savin Classification (women) to describe the degree of hair loss. Below are the pictures of the scales. The shaded areas show the type of hair loss that can be treated with the Laser Comb.

Norwood Hamilton
Classification – Males

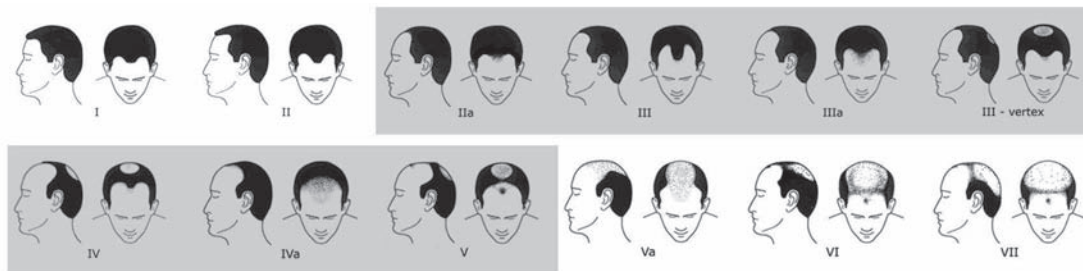

Ludwig-Savin
Classification – Females

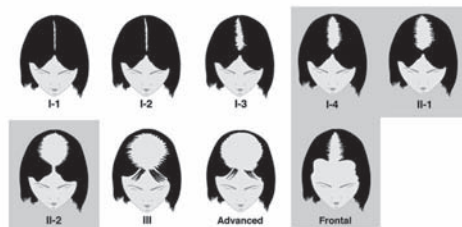

The Laser Comb has not been clinically tested for hair loss outside of the shaded areas.

The Laser Comb was only studied in light to medium skin tones (Types 1, 2, 3, 4) as described in the Fitzpatrick Scale below. If your hair loss and skin tone fall into these categories, the Laser Comb should work for you.

Fitzpatrick Skin Type Classification Guide:

A way to describe your skin type based on your complexion and what happens to your skin in sunlight.

Type 1	Clinically Tested to be effective for these skin types	Very white or freckled skin, always burns with sun exposure (very fair; often in people with red or blond hair and blue eyes)
Type 2		White skin, usually burns with sun exposure (fair; often in people with red or blond hair and blue, green, or hazel eyes)
Type 3		White or olive skin tone, sometimes burns with sun exposure (fair; seen in people with any hair or eye color)
Type 4		Brown skin, rarely burns with sun exposure (common in people of Mediterranean descent)
Type 5	Not Clinically Tested for these skin types	Dark brown skin, very rarely burns with sun exposure (common in people of Middle-Eastern descent)
Type 6		Black skin, rarely burns with sun exposure

Adapted from Fitzpatrick's Dermatology in General Medicine. McGraw-Hill Professional; 5th edition.

Please note that the Laser Comb was only tested on light to medium skin tones. This does not indicate that it will not work on darker skin tones.

WARNINGS AND CAUTIONS

WARNINGS

- ⚠ Do not stare directly at the red laser light or its reflection in a mirror because it can temporarily irritate your eyes. If you stare at the light for too long, it could harm your eyes. Also, never use a magnifying glass when using the Laser Comb because it can cause temporary irritation to your eyes.
- ⚠ Keep the Laser Comb out of the reach of children.
- ⚠ Do not drop the Laser Comb in water because you may get an electrical shock.
- ⚠ Do not use the Laser Comb near any heated surfaces because this could cause the device not to work properly and may cause an electrical shock.
- ⚠ Laser light - avoid direct eye exposure.

CAUTIONS

- ⚠ If you are taking any medications or products that make you sensitive to light, you should test the Laser Comb to make sure you do not have a reaction. Some of the things that can cause increased sensitivity to the Laser Comb could be medicines for colds, allergies, pain medicine and drugs to treat infections. **To test the Laser Comb for sensitivity**, turn the device on, hold the device on one spot, shining the light on your forearm. After two minutes, turn the device off and set aside. Look at the spot where you held the laser light over the next five minutes. If you see any reaction such as redness or rash, do not use the Laser Comb.
- ⚠ Only use the Laser Comb as directed. Do not try to change the way the Laser Comb operates. Changes may result in harmful laser exposure.
- ⚠ Never use acetone or other solvents on any part of the Laser Comb. Use of solvents on the window can block the laser light from coming through.

NOTICE: If any malfunction occurs, do not attempt to repair the Laser Comb yourself. This appliance has no user-serviceable parts. Use the Laser Comb only for its intended use as described in this manual. Contact Customer Service at info@anagenresearch.com.

RISKS

The FOLIGAIN® Laser Comb has been clinically tested on both men and women. In the clinical study for men, 5% said they had itching or tingling from using the Laser Comb. None of the women in the clinical study reported any of these side effects from using the Laser Comb. If you have itching, tingling or any other side effect that bothers you, stop using the Laser Comb.

Testing the Laser Comb for Sensitivity

Before beginning treatment, be sure to test the Laser Comb for sensitivity. Turn the device on and shine the light on your forearm, holding it in one spot. After two minutes, turn the device off and set aside. Look at the spot where you held the laser light over the next five minutes. If you see any reaction such as redness or rash, do not use the Laser Comb. If you don't experience any reaction, proceed with the treatment instructions on page 12 for the Advanced Laser Comb model or page 19 for the Professional Laser Comb model.

How to Monitor Device Activity

Your FOLIGAIN® Laser Comb is inspected and tested during production. When working properly, both the green and the red light indicators on the power and laser buttons turn on. A red laser light shines through the window at the bottom of the device. Do NOT look directly at the laser light. Shine the laser light onto the palm of your hand to see if the red laser light appears. If one or both of the light indicators is not lit, or if the laser light does not appear, contact the warranty provider for repair or replacement.

Device Failure Mode

The Laser Comb is designed to last for 10-15 years of continuous use. The laser light does not fade or diminish over time. At the end of the laser's life or if it fails for any other reason, the red laser light will no longer shine out of the device.

Expected Life of the FOLIGAIN® Laser Comb

The Laser Comb is a medical device manufactured under strict quality controls and validated through GMP (Good Manufacturing Practices). The Laser Comb contains high quality laser diodes designed to last at least 15,000 hours. Based on that rating and the durability of the components being used, the expected life of the Laser Comb is 10-15 years of regular use. Please take time to read this entire instruction manual before operating the Laser Comb.

Clinical studies were conducted at multiple locations in the United States. Participants in the study were males diagnosed with Androgenetic Alopecia, which includes hereditary hair loss, and had light to medium skin tones. Subjects were divided into two groups, Laser Comb Group and a Control Group. The Control Device looked and sounded like the Laser Comb device but did not include a laser light. The clinical studies adhered to all GCP (Good Clinical Practice) guidelines, were approved by an IRB (Institutional Review Board) and listed on www.clinicaltrials.gov.

During the 6-month clinical study, subjects treated their hair once per day, 3 times per week on non-consecutive days. The number of thick, normal, healthy hairs in the target zone was counted at week 1, week 16 and week 26.

Results at the 26-week visit showed that over 93% of the men saw an increase in hair counts (based on a minimum of 32 new hairs per square inch (6.45 square cm) being observed at the follow up visit).


No subjects in the study experienced any serious adverse events. In fact, the number and types of adverse events were similar in the Laser Comb and Control Groups.

Hair Count Change

In the clinical study, after 26 weeks of treatment with the Laser Comb devices versus the Control Device, the following increases in hair count were found.



**Mean Terminal Hair Count Changes
From Baseline in MALES
26 Weeks**

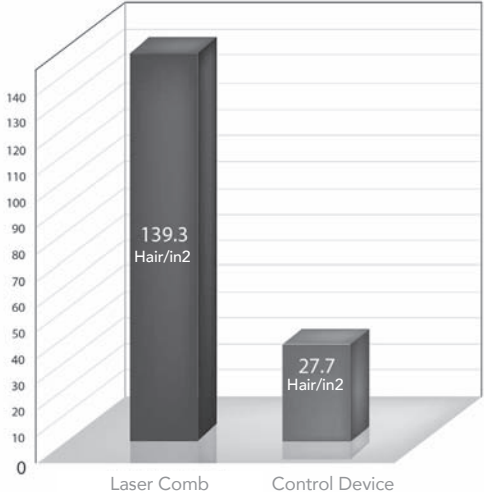

SUMMARY OF CLINICAL STUDY IN WOMEN

A clinical study was conducted at multiple locations in the United States. All participants in the study were females, diagnosed with Androgenetic Alopecia, which includes hereditary hair loss, and had light to medium skin tones. Subjects were divided into two groups. One group used the Laser Comb and the other group used a Control Device. The Control Device looked and sounded like the Laser Comb, but did not include a laser light. The clinical studies adhered to all GCP (Good Clinical Practice) guidelines, were approved by an IRB (Institutional Review Board) and listed on www.clinicaltrials.gov.

During the 6-month study, subjects treated their hair once per day, 3 times per week, on non-consecutive days. The number of thick, normal, healthy hairs in the target area was counted at Week 1, Week 16 and Week 26.


Results at 26 weeks showed that over 90% of the women saw increases in hair counts (based on a minimum of 32 new hairs per square inch (6.45 square cm) being observed at the follow up visit).

No subjects in the study experienced any serious adverse events. In fact, the number and types of adverse events were similar in the Laser Comb and control groups.


Hair Count Change

In the clinical study, after 26 weeks of treatment with the Laser Comb and a Control Device, the following increases in hair count were found.

**Mean Terminal Hair Count Changes
From Baseline in FEMALES
26 Weeks**

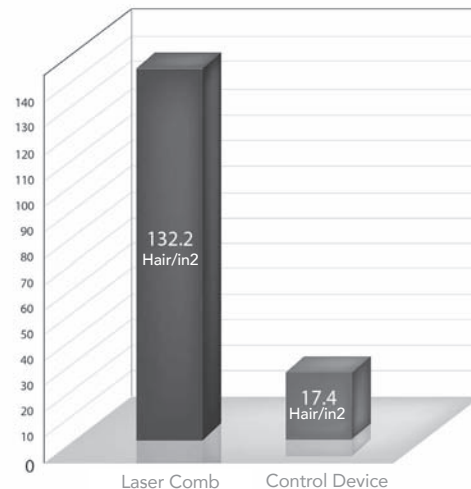

- A** POWER/LASER BUTTON
- B** BATTERY INDICATOR
- C** SAFETY LABEL
- D** TEETH
- E** BATTERY COVER TAB I
- F** LASER WINDOW
- G** WARNING LABEL
- H** CRADLE/STAND

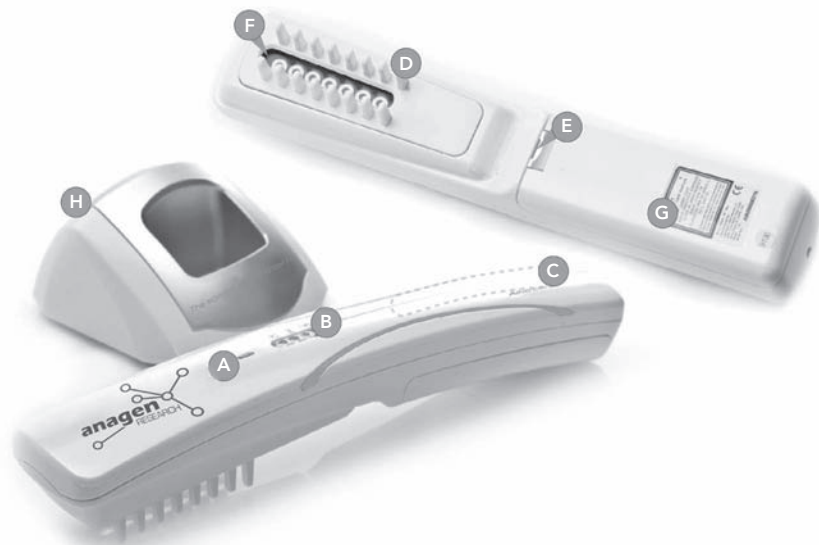


⚠ THIS SECTION APPLIES TO FOLIGAIN.L7x® ADVANCED LASER COMB MODEL ONLY


For FOLIGAIN.L12x® Professional Laser Combs refer to pages 16-23 for instructions

Your **Advanced Laser Comb** comes with four AA batteries already installed. To use the Laser Comb, you must remove the cardboard protective strip located in the battery compartment. Follow these steps:

1. **Verify the model of your device located on the handle.** If you purchased a Professional Laser Comb, refer to pages 16-23 of the User Manual.
2. Remove the blue label on the bottom of the Laser Comb.
3. Using a thin coin or your finger nail, depress the Battery Cover Tab (see Figure A) and remove battery door.
4. Open the battery compartment.

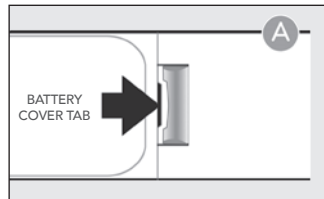

5. Locate and remove the cardboard protective strip (Figure B).
6. Replace Battery Cover by aligning bottom of Cover in grooves and then pressing the top of Cover until it snaps into place.
7. Your Laser Comb is now ready for use.

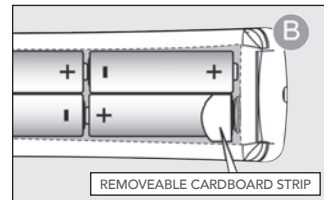

Replacing Batteries


The Laser Comb is designed to provide the optimum laser output even as the batteries age. When the batteries weaken to the point where they will no longer power the lasers to their optimum power level, the green battery indicator will flash from green to red and the lasers will flash.

This is an indication that the batteries need to be replaced. Use the following steps to replace the batteries.

Replacing batteries – continued

1. Using a thin coin or your finger nail, depress the Battery Cover Tab (see Figure A) and remove battery door.
2. Insert 4 AA batteries. Be sure to align batteries as shown on bottom of battery compartment (Figure C).
3. Replace Battery Cover by aligning bottom of Cover in grooves and then pressing the top of Cover until it snaps into place.

Do not use the Laser Comb if the device is damaged or defective; instead, contact HairMax Customer Service at info@hairmax.com

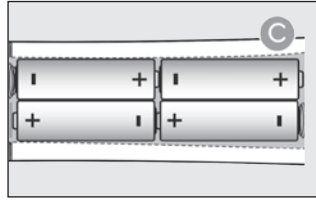

Proper Use

Read and follow all directions prior to use.


Use the Laser Comb for 15 minutes per treatment, three treatments per week.

The treatments should be spaced every other day, for example, on Monday, Wednesday and Friday. Do not use the Laser Comb more than 15 minutes per treatment day. Most users see benefits in as little as 16 weeks (based on clinical studies).



Slide the power/laser button to the middle or top position (Figure D).

Middle position activates the laser and vibrate feature. **Top position** is laser only. Shine the Laser Comb onto the palm of your hand. You will see that the laser light is on. (Do not look at the laser light directly). The Laser Comb is now ready for use.

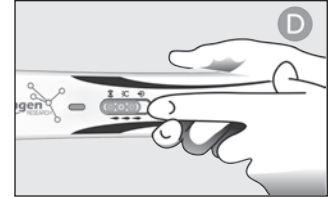

Treatment Instructions

Before beginning treatment, be sure to test the FOLIGAIN® Laser Comb for sensitivity (Page 5).



After you have turned on the Laser Comb, place it flat on your scalp so that both rows of teeth are touching your scalp (Figure E).




Start at your hairline. Move the Laser Comb from front to back of scalp in a slow combing motion, stopping every 1/2 inch (1cm) and pausing for 4 seconds (Figure F). Do not lift the Laser Comb off your scalp while combing from front to back. Next, lift the Laser Comb off your scalp and start at a new position at the hairline. Continue in this manner from front to back, then side to side, until you cover your entire scalp with multiple passes.

IMPORTANT: Do not lift the Laser Comb off your scalp while combing from front to back or side to side unless you are starting at a new position at the hairline. The Laser Comb features a pair of comb teeth that part the hair so the light can more effectively penetrate to the scalp. You will make a few complete passes over your scalp within recommended treatment time.

Repeat these passes within the required treatment of 15 minutes. If you have long hair, use your fingers to lift your hair out of the way allowing the Laser Comb to deliver the laser light to the scalp (Figure G).




Turn off the Laser Comb by sliding the power/laser button to the lowest position (Figure H). Store the device in the cradle/stand for future use.

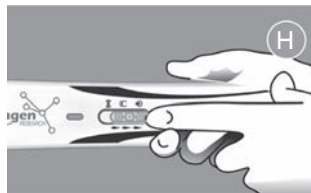

WHAT IS INCLUDED – FOLIGAIN.L7x® ADVANCED LASER COMB

Please Note:

- It is important to move the Laser Comb very slowly. Remember, you are bathing the scalp with laser light. If you move it too fast, the treatment will not achieve the optimal results to promote hair growth.
- Total time required for treatment with the Advanced Laser Comb is 15 minutes.

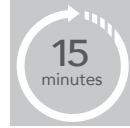

The Laser Comb is most effective when used correctly. In order to help time the movements, the Laser Comb will automatically vibrate every 4 seconds. **The vibration can be turned off** by sliding the power/laser button to the uppermost top position. The Laser Comb should be used for 15 minutes. Use a watch or clock to time your treatment.

What is included:

					
FOLIGAIN.L7x® Advanced Laser Comb	Retail Box	4 AA batteries	Cradle/Stand	User Manual	Quick Start Guide

The Laser Comb has advanced electronics to ensure that the device is working like it is supposed to and that the proper amount of laser energy is being delivered. If you are experiencing difficulties with your device, please try one of the tips below:

Your Laser Comb is turning off during use

- Make sure your Laser Comb has fresh batteries

Your Laser Comb will not turn off

- Slide the power switch to the lowest OFF position.

Your Laser Comb vibrates irregularly/doesn't vibrate

- Make sure the Power Switch is in the middle position.

Your Laser Comb vibrates, but you do not see the laser

- Make sure your Laser Comb has fresh batteries

You have trouble understanding any of these instructions

- Contact Customer Service at info@anagenresearch.com

SAFETY LABEL & LASER PARAMETERS	
Wavelength:	655nm +/- 10nm visible red light
Laser Power for classification:	<5mW CW
Divergence:	57 mrad
Number of lasers:	7
This Laser Product is designated as Class 3R during all procedures of operation	
Note:	The Laser Comb has no user-serviceable parts.
Location:	On the underside of the handle Model version varies



FOLIGAIN.L12x® PROFESSIONAL LASER COMB PARTS

- A** POWER BUTTON
- B** LASER BUTTON
- C** BATTERY INDICATOR
- D** INFORMATION DISPLAY
- E** POWERCORD
- F** SAFETY LABEL
- G** TEETH
- H** LASER WINDOW
- I** WARNING LABEL
- J** CHARGING CRADLE/STAND
- K** CHARGING CRADLE INDICATOR



⚠ THIS SECTION APPLIES TO FOLIGAIN.L12x® PROFESSIONAL LASER COMB MODELS ONLY

For FOLIGAIN.L7x® Advanced Laser Comb refer to pages 9-15 for instructions

Use the Laser Comb for 8 minutes per treatment, 3 treatments per week. The treatments should be spaced every other day, for example, on Monday, Wednesday and Friday. Do not use the Laser Comb more than 8 minutes per treatment day. Most users see benefits in as little as 16 weeks (based on clinical studies).




When you receive the Laser Comb, take all the parts out of the case. Remove the red safety label. The power button and the laser button are located under the red safety label. You must charge your Laser Comb for 2-3 hours before initial use.



Charging Options

Charging can be accomplished by one of the two following options:

OPTION 1: Charge the Laser Comb using the charging cradle (stand) and power cord. Prior to placing the Laser Comb in the cradle, make sure your Laser Comb is turned off and the teeth face away from you. Plug the power cord into the cradle and the other end into an electrical wall outlet. While the device is charging, you will see a red light on the charging cradle. After about 2-3 hours the light on the charging cradle will turn green, which means that the Laser Comb is fully charged and ready to use. During this charging option the light on the device will not light up.

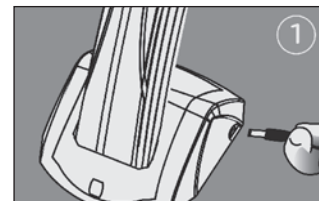

OPTION 2: Charge the Laser Comb by placing the power cord directly into the handle end of the Laser Comb and the other end into an electrical wall outlet. The battery indicator light on the top of the Laser Comb will show yellow when charging and turn green when fully charged.

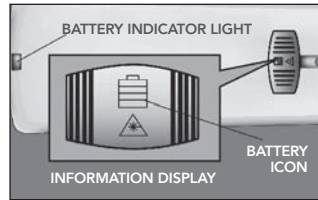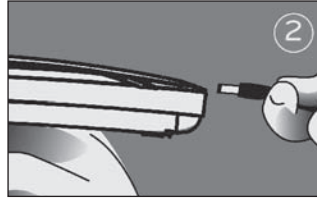


When the Laser Comb requires charging, the battery indicator light will flash red and the battery icon on the information display will appear. Charge your Laser Comb with charging option 1 or 2.

Proper Use

The Professional Laser Comb should be used for 8 minutes. Use a watch or clock to time your treatment.

First, press the power button and hold until power button light turns green. Second, press the laser button, which will light up red. That means the Laser Comb is ready for use.

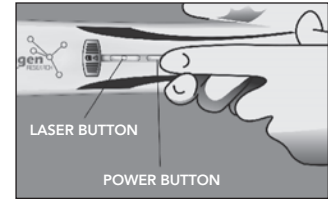


Shine the Laser Comb onto the palm of your hand and you will see that the laser light is on (do not look at the laser light directly). Do not use the Laser Comb if the device is damaged or defective; instead, contact HairMax Customer Service at info@hairmax.com.

Please Note:

If your device needs to be recharged during your treatment, simply insert the power cord directly into the end of the device and other end into an electrical wall outlet and continue with your treatment. When finished, allow the device to fully charge.

The information display on the top of the Laser Comb will show you 3 pieces of information:

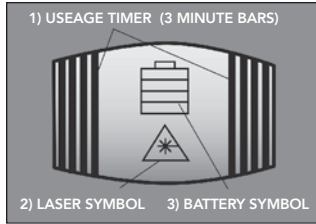


1. Usage Timer
2. Laser Symbol
3. Battery Symbol

- 1. Usage Timer** - Measures duration of treatment session.
This ensures all areas of your scalp are properly covered.
- 2. Laser Symbol** - To let you know that the laser is turned on.
- 3. Battery Symbol** - The battery icon will flash when a charge is needed.


Treatment Instructions

Before beginning treatment, be sure to test the FOLIGAIN® Laser Comb for sensitivity (Page 5).

After you have turned on the Laser Comb, place it flat on your scalp so that both rows of the teeth are touching your scalp (Figure A).

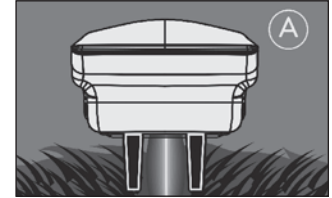

Start at your hairline. Move Laser Comb from front to back of scalp in a slow combing motion, pausing every 1/2 inch (1 cm) for 4 seconds (Figure B). Do not lift Laser Comb off your scalp while combing from front to back.

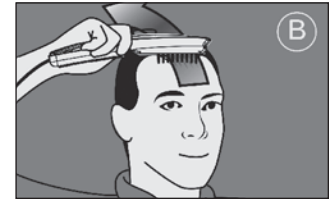

Next, lift the Laser Comb off your scalp and start at a new position at the hairline. Continue in this manner from front to back, then side to side, until you cover your entire scalp with multiple passes (Figure C).




Repeat these passes for the total treatment of 8 minutes. When the Laser Comb is on, it will beep and vibrate every 4 seconds. If you have long hair, use your fingers to lift your hair out of the way while using the Laser Comb to be sure that the laser energy reaches your scalp (Figure D).




IMPORTANT: Do not lift the Laser Comb off your head when moving every 4 seconds (unless you are checking the usage timer on the information display or until you complete the pass). The Laser Comb features a pair of comb teeth that part the hair so the light can more effectively reach the scalp.

When the Laser Comb is on, it will beep and vibrate every 4 seconds. This means that it's time to move the Laser Comb to the next area of the scalp. When you do this, do not lift the device off of your scalp.

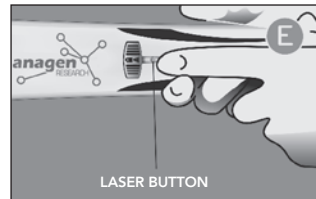

If you wish to **turn off the beep/ vibration** during the treatment, press the top laser button once (Figure E). Your Laser Comb and laser will remain on. To reactivate beep and vibration, simply press the laser button twice.

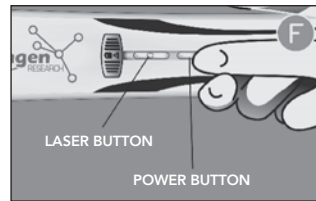

Turn off the Laser Comb by pressing and holding the power button (lower button) until green indicator light turns off (Figure F). Store the device in the charging cradle for future use. You may leave the power cord plugged into the charging cradle at all times.

Please Note:

- It is important to move the Laser Comb very slowly. Remember, you are bathing the scalp with laser light. If you move it too fast, the treatment will not achieve the optimal results to promote hair growth.
- Total time required for treatment with the Professional Laser Comb is 8 minutes.

What is included:

					
FOLIGAIN.L12x® Professional Laser Comb	Retail Box	Charging Cradle	Universal Adaptor	User Manual	Quick Start Guide

ANY ATTEMPT TO SELF-SERVICE YOUR DEVICE UNLESS INSTRUCTED BY A HAIRMAX REP. WILL VOID YOUR WARRANTY.

The Laser Comb has advanced electronics to ensure that the device is working like it is supposed to and that the proper amount of laser energy is being delivered. If you are experiencing difficulties with your device, please try one of the tips below:

Your Laser Comb will not turn on:

- Make sure your Laser Comb is fully charged.
- Insert the power cord into the Laser Comb, press the power button for 4 seconds and then press the laser button.
- Try another electrical wall outlet.

If your Laser Comb is turning off during use:

- Make sure your Laser Comb is fully charged.
- Insert the power cord into the Laser Comb, press the power button for 4 seconds and then press the laser button.
- Try another electrical wall outlet.

If your Laser Comb will not turn off:

- Press and hold power button for 4 seconds.
- If using power cord, unplug the Laser Comb and press power button for 4 seconds.

If your Laser Comb beeps/vibrates irregularly or doesn't beep/vibrate:

- If using power cord, unplug the Laser Comb. Plug in power cord again and press the power button, then the laser button.

If your Laser Comb beeps and vibrates, but you do not see the laser light:

- If using power cord, unplug Laser Comb. Plug in power cord again and press the power button, then the laser button.
- If you have trouble understanding any of these instructions Contact Customer Service at info@anagenresearch.com.

SAFETY LABEL & LASER PARAMETERS	
Wavelength:	655nm +/- 10nm visible red light
Laser Power for classification:	<5mW CW
Divergence:	57 mrad
Number of lasers:	12
This Laser Product is designated as Class 3R during all procedures of operation	
Note:	The Laser Comb has no user-serviceable parts.
Location:	On the underside of the handle Model version varies

Power Supply Specifications

Output = 6VDC 1.0A. Plug: 9mm x 3.5mm, center positive.

The power supply for the Laser Comb can be used internationally. This power supply will adapt to all power supplies from 110 volt/60 cycle (US) to 240 volt/50 cycle. An adapter specific to your country is included with your Laser Comb.




CLEANING YOUR FOLIGAIN® LASER COMB

- Keep your Laser Comb clean.
- Turn off the Laser Comb by sliding the power button to the lowest position (FOLIGAIN.L7x® Advanced) or by pressing and holding the power button (lower button) until green light turns off (FOLIGAIN.L12x® Professional).
- Unplug the power cord from the device and from the electrical wall outlet.
- Clean the laser window with a soft, slightly moistened cloth.
- To clean the comb teeth, gently remove any hairs or dust with a slightly moistened towel.
- Don't use any acetone, soap or detergents. Use warm water only.

- To clean the teeth more thoroughly, gently pull out the comb teeth from the center. Then you can separately wash them using soap and warm water.

- After you dry the comb teeth, put them back on the device by

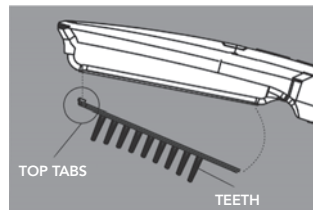

inserting the top tabs into place and press the bottom of the teeth with your thumb until they click into place. If you need to replace the teeth, please contact info@anagenresearch.com.

Electrical Safety

Basic safety precautions should be observed when using this, or any other electrical device, especially when children are present.

⚠ WARNING: Keep away from water. Contact with water may result in an electrical shock.

To reduce the risk of hazard by electric shock:

- Always unplug the Laser Comb from the power outlet immediately after using.
- Do not use while bathing.
- Do not place or store the Laser Comb where it can fall or be pulled into a tub or sink. Do not expose the Laser Comb to running water, or drop into water or any other liquid. If the Laser Comb falls into water, unplug it immediately from the power outlet.
- If using the Laser Comb near water, the use of an electrical outlet protected by a Ground Fault Circuit Interrupter is recommended.
- Do not dispose of the Laser Comb in a fire or store in high temperatures greater than 140°F (60°C).

WARRANTY

ANY ATTEMPT TO SELF-SERVICE YOUR DEVICE UNLESS INSTRUCTED BY A HAIRMAX REP. WILL VOID YOUR WARRANTY.

THE FOLIGAIN® LASER COMB POWERED BY HAIRMAX IS WARRANTED BY LEXINGTON RESEARCH, LLC ("MANUFACTURER") AGAINST DEFECTS IN MATERIALS OR WORKMANSHIP FOR 24 MONTHS FROM DATE OF PURCHASE. IF THIS PRODUCT IS DETERMINED TO BE DEFECTIVE, LEXINGTON WILL REPAIR OR REPLACE IT, AT ITS OPTION, AT NO CHARGE DURING THE WARRANTY PERIOD.

AFTER THE WARRANTY PERIOD, THE USER MUST PAY FOR ALL PARTS AND LABOR CHARGES. IF THE PRODUCT IS REPAIRED OR REPLACED, EITHER WITHIN OR OUTSIDE THE WARRANTY PERIOD, LEXINGTON RESERVES THE RIGHT TO PROVIDE EQUIVALENT PARTS IN DIFFERENT COLORS TO YOUR ORIGINAL PRODUCT.

THIS WARRANTY DOES NOT COVER DAMAGES CAUSED BY ACTS OF GOD, MISUSE, NEGLIGENCE, ACCIDENT, MODIFICATION OF, OR TO ANY PART OF THE PRODUCT. THIS WARRANTY DOES NOT COVER DAMAGE DUE TO IMPROPER OPERATION OR MAINTENANCE, CONNECTION TO IMPROPER VOLTAGE SUPPLY, OR ATTEMPTED REPAIR BY ANYONE OTHER THAN LEXINGTON. THIS WARRANTY IS VOID IF PRODUCT IS TAMPERED WITH.

TO OBTAIN WARRANTY SERVICE, YOU MUST DELIVER THE PRODUCT FREIGHT PREPAID, IN EITHER ITS ORIGINAL PACKAGING OR PACKAGING AFFORDING AN EQUAL DEGREE OF PROTECTION, TOGETHER WITH YOUR NAME AND ADDRESS, TO:

Lexington Research, LLC, 777 Yamato Road, Suite 105, Boca Raton, FL 33431, USA

LIMITATION OF LIABILITY; INDEMNITY: EXCEPT AS EXPRESSLY SET FORTH ABOVE, LEXINGTON RESEARCH, LLC SHALL NOT BE LIABLE FOR ANY SPECIAL, INDIRECT, OR CONSEQUENTIAL DAMAGES THAT ARE DIRECTLY OR INDIRECTLY RELATED TO THE USE OF, OR THE INABILITY TO USE, THE FOLIGAIN® LASER COMB POWERED BY HAIRMAX, EVEN IF LEXINGTON RESEARCH, LLC HAS BEEN ADVISED OF THE POSSIBILITY OF SUCH DAMAGES WHETHER IN AN ACTION OF CONTRACT, NEGLIGENCE OR OTHER TORTUOUS ACTION. SOME STATES DO NOT ALLOW THE EXCLUSION OR LIMITATION OF INCIDENTAL OR CONSEQUENTIAL DAMAGES, SO THE ABOVE LIMITATION OR EXCLUSION MAY NOT APPLY TO YOU.

GOVERNING LAW: THE USAGE OF THE FOLIGAIN® LASER COMB AND THE INTERPRETATION OF ITS TERMS SHALL BE GOVERNED BY AND CONSTRUED IN ACCORDANCE WITH THE LAWS OF THE STATE OF FLORIDA APPLICABLE TO AGREEMENTS NEGOTIATED, EXECUTED AND PERFORMED WHOLLY AND SOLELY WITHIN THE STATE OF FLORIDA. THE PARTIES AGREE THAT THE STATE OF FLORIDA'S CONFLICTS OF LAWS RULES SHALL NOT APPLY TO THIS AGREEMENT OR ANY TRANSACTION PERFORMED BY OR UNDER THIS AGREEMENT.

JURISDICTION AND VENUE: THE PARTIES IRREVOCABLY SUBMIT AND CONSENT TO THE SOLE AND EXCLUSIVE JURISDICTION AND VENUE OF THE FLORIDA STATE COURTS IN AND FOR PALM BEACH COUNTY, FLORIDA, U.S.A. AND THE FEDERAL COURTS IN AND FOR THE SOUTHERN DISTRICT OF FLORIDA, U.S.A. THE PARTIES AGREE NOT TO RAISE THE DEFENSE OF FORUM NON CONVENIENS.

AUTHORIZED REPRESENTATIVE IN EUROPE (REGULATORY AFFAIRS ONLY):

EMERGO EUROPE, MOLENSTRAAT 15, 2513 BH THE HAGUE, THE NETHERLANDS. TEL: (31) (0) 70 3t5-8570 FAX: (31) (0) 70-345-7299

To validate your warranty please register online at www.hairmax.com/warranty.

WARRANTY CARD

The FOLIGAIN® Laser Comb has a 2 year warranty against defects in materials or workmanship from Lexington Research, LLC. Items that are damaged or abused are not covered under the warranty and will be charged at the regular repair rates. **ANY ATTEMPT TO SELF-SERVICE YOUR DEVICE UNLESS INSTRUCTED BY A HAIRMAX REP. WILL VOID YOUR WARRANTY.** Anagen Research®, LLC. may extend this warranty period under certain product failures. Please refer to the User Manual for complete warranty terms. We guarantee your satisfaction.

To validate your warranty card register online at www.hairmax.com/warranty

or complete and mail this warranty card.

Name: _____ Serial Number: _____

*located on yellow safety label,
on the bottom of the device*

Date of Purchase: _____

Address: _____

City: _____ State: _____ ZIP/Postal Code: _____

Country: _____ Email: _____

Where did you purchase this product (Reseller Name)? _____

What factors influenced your purchase of this product? _____

**Please fill out this form, detach it, place it in an envelope, or fold out and mail it directly to us with a postage stamp.*



Staple or
tape here

Return Address:

Place
stamp
here

Lexington Research, LLC
777 Yamato Road, Suite 105
Boca Raton, FL 33431





Anagen Research® Customer Service
Email: info@anagenresearch.com
www.anagenresearch.com

Mailing Address:
8560 West Sunset Boulevard, Suite 500
West Hollywood, CA 90069 USA





FOR INTERNATIONAL TRANSLATIONS VISIT:

www.anagenresearch.com

- FR Pour une traduction internationale, visitez: www.anagenresearch.com
- ES Para traducciones internacionales, por favor visite: www.anagenresearch.com
- IT Per traduzioni internazionali, per favor visitate: www.anagenresearch.com
- DE Für internationale Übersetzungen besuchen Sie bitte: www.anagenresearch.com
- PT Para as traduções internacionais, visite: www.anagenresearch.com