



Beautistyle International Corp.

3F, No. 10, Aly. 6, Ln. 45, Baoxing Rd. Xindian City, Taipei County, Taiwan, R.O.C.

TEL: +886-2-2918-9958

FAX: +886-2-2918-2068

www.beautistyle.com.tw

All rights reserved © Beautistyle International Corp.

Version 1, Dec 2009/11/28

Ozone Beauty Kit

AR-389

User Manual



Protect Your Beauty
Oxygenate the Skin

Content

Introduction / Product Specifications	p2
Product Diagram/ Functions.....	p3
Home Beauty/ Home Health	p4
Restrictions /Precautions	p5
Warranty Card.....	p6

§ Revitalize your skin with the Ozone Beauty Kit

The Ozone Beauty Kit uses the "Plate-type Silent Electric Discharge" method of super-oxygenation to convert oxygen in the air from $3O_2$ to $2O_3$. The resultant ozone is a clean and refreshing gas that has anti-bacterial, detoxification, de-odorizing and bleaching effects. It is able to break down grime and balance oils. The blue light cleanses, whitens and relaxes so the facial skin no longer looks greasy and resists acne. You will rediscover your healthy and beautiful looking skin in no time. (*extract:Ozone Therapy – A New Total Approach to Health*)

§ Product Specifications

Product Name: Ozone Beauty Kit

Weight: 136mg


Dimensions: 19.5 x 6.5 x 5cm

Power Supply: 100 ~ 240V

§ Parts



Protective Cover

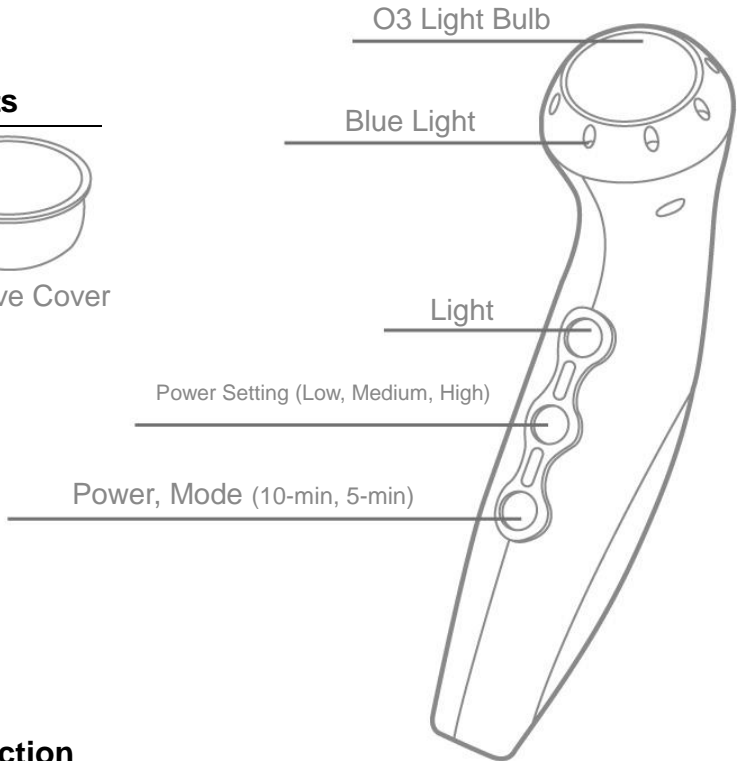

§ Function

1. Plug the Ozone Beauty Kit into a power supply and press the POWER button to turn it on.
2. Mode Selection: The Power key is also used to set the Mode. Press once for the green light to come on and the kit will run for 10 minutes. Press again for the red light and the kit will run for 15 minutes. Press again to turn off.
3. Strength Selection: When the kit comes on the default strength is Low (green light). Press the Setting button for Medium (yellow light) and again for High (red light). Press again to revert back to Low (green light).
4. Light: The Light button turns the Blue Light on and off.
5. When the Ozone Beauty Kit is on and in contact with the skin, it emits a sizzling sound, ultra-violet light (only visible in the dark) and a grassy scent. Do not be concerned as this is normal.
6. The Blue Light is for localized revitalization after using the Ozone Beauty Kit.

§ Home Beauty

Step 1 – Stable Ozone

When using this product, turn on for 30 seconds to let the electrodes break down oxygen in the air from $3O_2$ into $2O_3$ and produce stable ozone.

Step 2 – Deep Cleanse


After cleansing the face, use this product at the T and V parts of the face to fully cleanse the epidermis, prevent aging, whiten the skin and promote metabolism.

Avoid touching acne/pimples on the face directly with this product. Run the kit lightly over the surrounding areas and wait for the affected part to scab over before treatment.



§ Home Care

For sensitive, itchy or allergic skin, as ozone has an anti-bacterial effect, touch the affected area lightly with the O_3 light bulb to sterilize the skin and prevent itchiness.

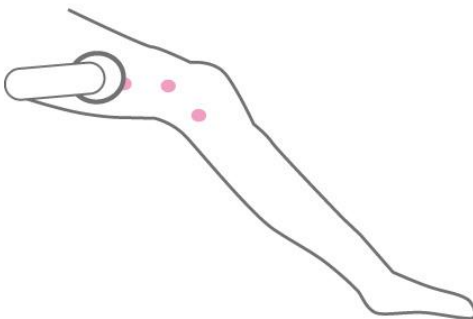

§ Restrictions

1. Do not use if you have a heart disease or a pacemaker fitted.
2. Do not use if you have a malignant tumor.
3. Do not use during menstruation or pregnancy.
4. Do not use if you have skin disease, sores or open wounds.
5. When using the blue light, please close your eyes or wear a blind-fold as the strong light can cause discomfort to the eyes.

§ Precautions

1. Avoid using the beauty kit over the eyes. Do not use continuously on the same location for more than 3 minutes.
2. A mild stinging sensation is normal during use. The strength setting can be adjusted depending on personal preferences.
3. Clean with a soft (cotton) cloth after use. Do not wash directly with water.
4. Store in a dry location. Do not leave kit in the bathroom.
5. The O₃ light bulb is made from a glass material and may shatter under strong impact.
6. Do not attempt to disassemble or modify the product. If the product malfunctions, please send it back to our company for technical support.


Warranty Card

Customer Details

Name: _____ Tel: _____

Address: _____

Distributor/Retailer:

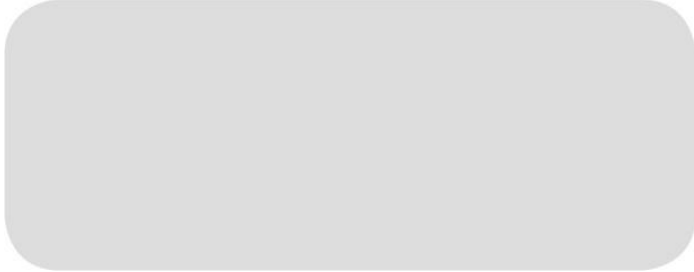

Date Purchased:

Please retain this user manual and warranty card to protect your rights as a consumer.

(1) This product warranty is for one year from the date of sale and free repairs are provided by Beautistyle under normal conditions of use. Beautistyle may charge for labor and parts under the following circumstances:

- i. Damage is due to improper use, operation or modification.
- ii. Damage due to factors beyond human control.
- iii. If product is out of warranty.

(2) Please keep this warranty card in a safe place as no replacement is available if lost.

(3) This warranty card is only valid if stamped with a company seal.

(4) Consumables are not covered by this warranty.

(5) Please keep the packaging and accessories clean and intact.

All parts included with this product were fully tested before sale. If your product is defective or damaged, please return it for replacement within 7 days of purchase. Replacements will not be possible if the product has any missing or damaged parts.