1-877-518-2323

ask@bmrbeautyseries.com

SDI Inc.,
A division of Bio-Medical Research Ltd.,
PO Box 5179
Hoboken, NJ 07030

Designed by & Manufactured for:
Bio-Medical Research Ltd., Parkmore Business Park West, Galway, Ireland.
Copyright: © 2011, Bio-Medical Research Ltd. All Rights Reserved
Part No: 2400-0330 Rev.: 1 Date of Issue: 2/11

*Arm
Lift*

INSTRUCTION MANUAL

**MODE D'EMPLOI
MANUAL DE USARIO**

www.bmrbeautyseries.com

CONTRA-INDICATIONS, WARNINGS & PRECAUTIONS

BMR ARM LIFT is indicated for the toning of the triceps muscles of the upper arms. Your product is suitable for use by all healthy adults; however, as with other forms of exercise, some care is needed when using it. Always follow the guidelines below and read the instruction manual before use. Some of the points below are gender-specific.

CONTRA-INDICATIONS:

Do not use if:

- You have an electronic implant (e.g. cardiac pacemaker or defibrillator - as your ACTIVATOR may interfere with the proper functioning of the implanted stimulator), or if you suffer from any other heart problem.
- You have any other form of medical implant. The safety has not been established for the use of electronic muscle stimulators over medical implants.
- You are pregnant. The safety of electronic muscle stimulation (EMS) over the pregnant uterus has not been established.
- You suffer from cancer. The effects of electronic stimulation on cancerous tissue are unknown.
- You are under medical supervision for cognitive dysfunction as you may not be able to comply with safety instructions.
- The ACTIVATOR is in close proximity (e.g. 1 yard) to shortwave or microwave therapy equipment or you are connected to high-frequency surgical equipment (e.g. diathermy), because of the risk of device malfunction due to electronic interference.
- Wearing the device necessitates placement over areas at which drugs/medicines are administered (short-term or long-term) by injection (e.g. hormone treatments).

WARNINGS:

- The long-term effects of chronic electronic stimulation are unknown.
- This product should not be used while sleeping, because you cannot react to safety information.

When applying the GelPads™ and garments, always remember to:

- Place the GelPads and garments ONLY as indicated in this manual.
- Not place the GelPads or garments on the front or sides of the neck, across or through the heart (i.e. one GelPad on the front of the chest and one on the back), in the genital region or on the head, because of the risk of stimulating inappropriate muscles and organs.
- Avoid application of GelPads near the thorax as it may increase the risk of cardiac fibrillation.
- Avoid any recent scars, broken or inflamed skin, areas of infection or susceptibility to acne, thrombosis or other vascular problems (e.g. varicose veins), or any parts of the body where feeling is limited.
- Avoid areas of injury or restricted movement (e.g. fractures or sprains).
- Avoid placing the GelPads over metal implants.
- BMR ARM LIFT IS not intended for medical use, for the treatment of any medical condition, or for any permanent physical changes.

PRECAUTIONS:

Please read these precautions before using BMR ARM LIFT:

- You should only use this product while in a relaxed, seated position, with your arms bent slightly and resting in your lap.
- Do not straighten your arms during the session, as this can lead to an ineffective stimulation and may cause some discomfort.
- Never use BMR ARM LIFT while involved in any activity which includes using your arms or hands.
- Your arms are likely to move involuntarily during a session.
- Keep the ACTIVATOR in your hands during each session.
- Do not fasten the straps too tightly around your arms.
- If the garment leads are damaged or exposed, they should be replaced immediately.
- Maximum use of a BMR ARM LIFT product is one session per day per body area.
- Caution should be used for patients with suspected or diagnosed epilepsy.
- Caution should be used in the following:
 - When there is a tendency to bleed internally following an injury or fracture.
 - Following recent surgical procedures when electronic stimulation may disrupt the healing process
 - Over areas of skin which lack normal sensation

- This product is not intended for medical use, for the treatment of any medical condition, or for any permanent physical changes.

Get your doctor's permission before using this product if:

- You have any serious illness or injury not mentioned in this guide.
- You have recently undergone a surgical procedure.
- You take insulin for diabetes.
- You want to use it on a young child.
- You suffer from muscle or joint problems.

Adverse reactions:

- Skin irritations and burns have been reported with the use of powered muscle stimulators.
- On very rare occasions, first-time users of EMS have reported feeling light-headed or faint. We recommend that you use the product while seated until you become accustomed to the sensation.

To reposition GelPads during a session:

- Always pause the program currently running, unfasten the strap and then refasten it around your arm after the GelPads' positions have been adjusted.

After strenuous exercise or exertion:

- Always use a lower intensity to avoid muscle fatigue.

Contact BMR or an authorized dealer if:

- Your BMR ARM LIFT is not working correctly. Do not use in the meantime.
- You experience any irritation, skin reaction, hypersensitivity, or other adverse reaction. You should, however, note that some reddening of the skin under the garments is quite normal during and for a short time after an exercise session.

Note:

- An effective session should not cause undue discomfort.

Important:

- For first time users, muscle stimulation can be an unusual sensation. We recommend that you begin in a seated position with low stimulation intensity settings to familiarize yourself with the sensation before progressing to higher intensity settings.
- Keep your BMR ARM LIFT out of the reach of children.
- The studs and GelPads must not be connected to other objects.
- Do not use your BMR ARM LIFT at the same time as any other device which transfers an electrical current into the body (e.g., another muscle stimulator).
- Do not touch the GelPads or metal studs while the unit is switched on.
- Do not use while driving, operating machinery or cycling.
- The product is for single person use only.
- Cease using BMR ARM LIFT if you are feeling light-headed or faint. Consult your doctor if this happens.
- Do not over exert yourself when using muscle stimulation. Any workout should be at a level comfortable for you.
- Although compliant with applicable EMC requirements, this device may still interfere with more sensitive equipment, please move away or turn off.
- BMR will not accept responsibility if the guidelines and instructions supplied with this product are not followed.

Note: If you are in any doubt about using this product for any reason, please consult your doctor before use.

WELCOME

Congratulations on purchasing your BMR ARM LIFT. Please read your user manual fully before using this product.

This revolutionary new concept allows a single hand-held ACTIVATOR (sold separately) to be attached to a range of garments which tone different areas of the body, saving you money on your next BMR LIFT purchase. BMR ARM LIFT uses clinically proven EMS technology to exercise the triceps muscles of the upper arms. A signal is sent to the muscles and as that signal strength increases, the muscles contract. The product should be used while relaxing in the comfort of your own home.

Use your BMR ARM LIFT 5 times a week for 10-20 mins and in just 4 weeks, you can expect to see more toned arms and have a renewed confidence in wearing sleeveless tops and dresses!

For best results we recommend that you use your BMR ARM LIFT in conjunction with a normal, healthy diet and exercise.

BMR ARM LIFT is an electronic muscle stimulator intended for the enhancement of the performance of normal muscle.

If you have any difficulty setting up your product or have any other questions or concerns about using your BMR ARM LIFT, contact the Careline at:

1-877-518-2323

Email: ask@bmrbeautyseries.com

www.bmrbeautyseries.com

CARING FOR YOUR GARMENTS

Your garments may be cleaned using a lightly dampened sponge, but you must first remove the ACTIVATOR and GelPads. Always follow the cleaning instructions when cleaning the garments.

Never machine-wash or hand-wash your garments. Always clean them using a lightly dampened cloth or sponge.

Do not use bleach when washing the garments.

Do not dry clean your garments.


Do not tumble dry your garments. Dry them on a flat surface. Do not dry them over anything hot. (e.g. a radiator) as they contain plastic parts. Ensure the garments are completely dry before using them again.

The garments should never be ironed.

INTRODUCING YOUR BMR ARM LIFT

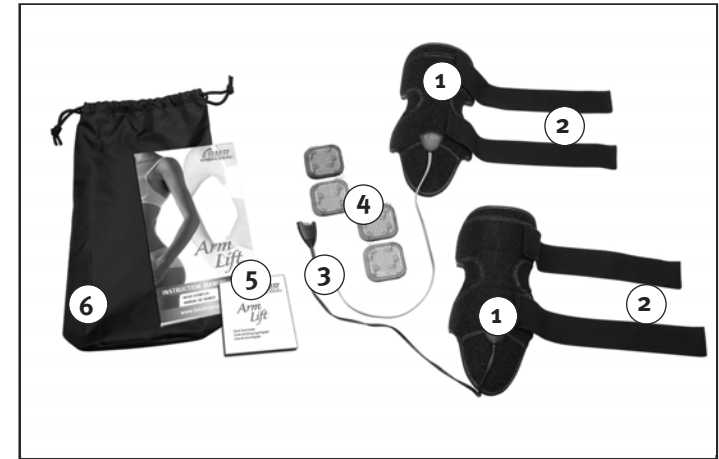You should find the following parts inside your pack.

1. BMR ARM LIFT garments (x 2):

The arm garments have been ergonomically designed and are fastened around your upper arm using the straps.

2. Straps (x 4):

These are used to hold the garments in place. They will fit an arm circumference range between 8.5" - 14.5".

3. Unit Connector:

This connects the garments to the ACTIVATOR (sold separately). Make sure the ACTIVATOR and garments are fully connected before use.

IMPORTANT: Do not disconnect the ACTIVATOR from the garments until it has been switched off.

4. Pack of four adhesive GelPads:

These GelPads are placed onto the garments before use. When you wrap the garments around your arms, these pads must be in contact with your skin so that signals can be sent to your muscles.

5. Instruction Literature:

This is a detailed guide to using your BMR ARM LIFT and an additional quick start guide to help get you started.

6. Carry Pouch:

Use this pouch to store your BMR ARM LIFT between toning sessions.

ACTIVATOR CONTROLS

1. On/Off Button (⏻)

Press and hold this button for 2 seconds to turn your ACTIVATOR on or off. You may also pause the exercise by briefly pressing this button.

2. Increase Toning Intensity (▲)

Press and hold these buttons to increase the intensity of the exercise. The right-hand button controls the garment on your right arm. The left-hand button controls the garment on your left arm.

3. Decrease Toning Intensity (▼)

Press and hold these buttons to decrease the intensity of the signal to either arm.

4. Program Button (P)

Press this button to select the toning program you wish to use. There are 3 programs, which are explained later in this manual (page 10).

5. Information Button (i)

Press this button to see your training information, such as the maximum intensity for your current or last session and the number of sessions you have completed.


ACTIVATOR DISPLAY

- Toning intensity from the left arm garment.
- Toning intensity from the right arm garment.
- Counts down the time left in the current session. Also displays error messages (see page 11 & 13).
- Appears when the exercise has been paused.
- Indicates that the sound is off.
- Keylock is active - you cannot change the intensity levels or the program accidentally.
- This symbol indicates a poor contact between the ACTIVATOR and garments or between the GelPads and the skin (see page 13 for more details).
- Shows the battery power remaining.
- Number of sessions completed.
- Shows the program currently running (1 - 3).
- This symbol appears during each contraction phase.
- Maximum intensity for the current/ last session.


HOW DOES MUSCLE STIMULATION TECHNOLOGY WORK?

BMR ARM LIFT uses electrical muscle stimulation (EMS) technology. This technology has been used in hospitals and by physical therapists for over 50 years. Here is a summary of how EMS technology works.

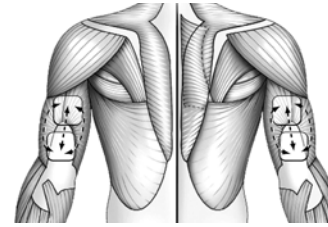Step 1

Signals are sent between the GelPads. The signals switch on the nerves which control the triceps muscles.

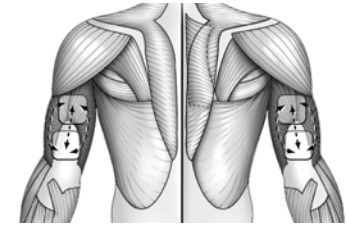
Step 2

As the signal strength increases, the nerves, in turn, cause your muscles to contract and relax rhythmically.

SETTING UP YOUR BMR ARM LIFT


Step 1 - Charge the battery. Insert the charger into a socket and then connect it to your ACTIVATOR (Fig. a). It will take approximately 2-3 hours to fully charge the battery. Your ACTIVATOR is fully charged when all three sections of the battery icon are full. This icon will flash when the battery is low and needs to be recharged.
NOTE: Do not leave your ACTIVATOR connected to the charger when the battery is fully charged.

Fig. a

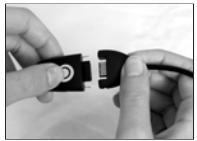
Step 2 - Place the GelPads on the garments. Place the garments so that the four silver studs are facing you. Remove the GelPads from their pack. You will notice that one side of each GelPad has a grid pattern on it while the other side is plain black.
NOTE: You should only use BMR ARM LIFT GelPads with this product

Fig. b

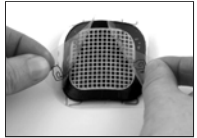(i) Remove the covers from the **patterned side** of the four GelPads (Fig. b).

Fig. c

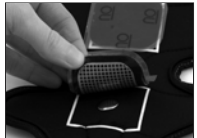
(ii) Place the **patterned side** of the GelPads onto the silver studs within the white squares (Fig. c).

IMPORTANT:

Before you progress, make sure that all of the silver studs are completely covered by the GelPads. Press the edges of each GelPad firmly onto the material before use.

Fig. d

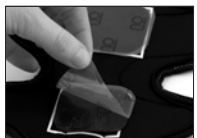Step 3 - Remove the covers from the black side of all four GelPads (Fig. d). Do not throw these covers away as you will need to put them back onto the GelPads at the end of your session.

Step 4 - Positioning the garments.

(i) Attach two straps to one side of each garment (fig. e).

(ii) Take the garment which is connected to the **GOLD** lead and place the point of your **RIGHT** elbow directly onto the connector symbol (🔌) on the inside of that garment. Then press the garment onto the back of your arm (Fig. f). Fasten that garment firmly in place using the straps (Fig. g).

(iii) Repeat step (ii), placing the garment connected to the **BLACK** lead onto your **LEFT** arm.

IMPORTANT:

You may need to adjust the position of the garments slightly until you find the most comfortable position. However, the positioning will become much easier after a few sessions.

Step 5 - Connect the ACTIVATOR to the garments. The **ACTIVATOR** and garments must be connected properly for your **BMR ARM LIFT** to work. 'Click' the **ACTIVATOR** connector into the garment connector as shown in Fig. h.

NOTE:

For the best results, we recommend that you complete five sessions per week. You should, however, only do one session per day, as this allows your muscles to recuperate.

Fig. e

Fig. f

Fig. g

Fig. h

Muscle toning can be an unusual sensation, but a pleasant one. It may tickle to begin with. This progresses to a smooth muscle contraction as the toning intensity increases. We recommend that you use the product while seated until you become accustomed to the sensation. To get started, follow these simple steps:

Step 1: Press and hold the on/off button (🔌) for two seconds to switch on your **ACTIVATOR** (Fig. a).

Step 2: Select the program you wish to use by pressing the program button (P - Fig. b). Then, to start the program, press and hold the increase intensity buttons (▲ - Fig. c) until you feel your muscles contracting. Always choose an intensity level at which you feel a strong but comfortable contraction. This will vary from person to person. The increasing intensity level is displayed on the screen.

Step 3: Continue increasing the intensity throughout the session if possible. You will feel the muscle contractions getting stronger as the intensity increases. Remember, the harder you work your muscles the better the results you will see, but always reduce the intensity if you feel any discomfort. Your unit's intensity range is 0-99.

PAUSING A SESSION

If you wish to pause a session before it is finished, simply press the on/off button briefly. The display will show that the session is paused (||). To resume the session, briefly press the on/off button again. The display will return to its normal mode.

Step 4: SESSION OVER

At the end of the training session your **BMR ARM LIFT** will stop automatically. However, you may switch off your **ACTIVATOR** at any time during a session by pressing and holding the on/off button for two seconds. You should see the display turn off.

Step 5: Remove the garments and put the covers back on the black side of the GelPads. Store everything in the pouch until your next session.

NOTE: Do not disconnect the **ACTIVATOR** from the garments without first switching it off, as this may result in an error on your display.

Well done on completing your first training session!

Fig. a

Fig. b

Fig. c

GETTING RESULTS WITH BMR ARM LIFT

- 1. Use the strongest intensity you can** - the harder you work your muscles, the better the results you will see, but remember, stimulation should never be uncomfortable.
- 2. Use your BMR ARM LIFT as part of a normal, healthy lifestyle.** Improve your diet and try to increase the amount of exercise you do. Other BMR LIFT garments are available to help tone and firm other areas of your body.
- 3. Your BMR ARM LIFT has 3 toning programs which are listed in the table below.** The ACTIVATOR automatically progresses from one program to the next after a set number of sessions. The display shows the program currently selected.

Program	Duration	Training Level
Beginner	10:00	Light
Intermediate	15:00	Moderate
Advanced	20:00	Strong

You cannot change a program during a session. First switch your ACTIVATOR off and on again. Then select your program using the program button.

- 4. To help you get the most from BMR ARM LIFT, you should follow the 30-day plan.** To do this, use BMR ARM LIFT 5 times a week for 4 weeks and record the highest intensity levels in your diary (the diary is inside the back cover). Each week you have two rest days. The diary below shows one person's plan, recorded during scientific testing over a 30-day period. The intensity level will vary from person to person.

Try to increase your toning intensity regularly. But remember, do not over-exert yourself. Any workout should be at a level comfortable for you.

30-DAY PLAN

Week	Day 1	Day 2	Day 3	Day 4	Day 5	Day 6	Day 7
1	(20 20)	rest	(25 25)	rest	(30 30)	(40 40)	(40 40)
2	rest	(50 50)	(55 55)	rest	(60 60)	(65 65)	(65 65)
3	rest	(65 65)	(70 70)	rest	(70 70)	(70 70)	(75 75)
4	rest	(75 75)	(75 75)	rest	(80 80)	(80 80)	(80 80)

Each box shows the highest toning intensity for either arm during each session.

Fill-in your diary every time you complete a toning session. Leave the diary somewhere you will see it regularly to motivate you to reach your goal.

ADDITIONAL ACTIVATOR FUNCTIONS

Information Button (i)

Press the information button any time during a session to see the highest intensity level you have reached for that session. Pressing the information button twice shows you the total number of sessions you have completed (Fig. a). While pressing the information button three times identifies the product type you are using.

Fig. a

Mute Function (M)

If you want to switch off your ACTIVATOR's sound effects, press and hold the program button for two seconds (Fig. b). The mute function remains active indefinitely unless manually changed. Deactivate the mute function by again pressing the program button for two seconds.

Fig. b

Keylock Function (🔒)

If you find a comfortable exercise intensity, press and hold the information button for two seconds to lock that intensity (Fig. c). This function only remains active during the session in which it is activated. If you wish to increase the intensity further, deactivate the keylock function by again pressing the information button for two seconds.

Fig. c

Error Messages

In the unlikely event of your ACTIVATOR developing a problem, "Err" will appear on your ACTIVATOR display (Fig. d). If this occurs you should switch the ACTIVATOR off and then switch it back on again. It should now operate properly. If the problem persists, please call the Careline for further assistance:

1-877-518-2323

Email: ask@bmrbeautyseries.com

www.bmrbeautyseries.com

Fig. d

Battery Power / Replacing the Battery

The battery icon on the display will flash when the battery is low and needs to be recharged. After a period of time, you may find that the battery will only power your ACTIVATOR for a limited number of sessions, requiring you to charge the battery more frequently. This is because rechargeable batteries have a limited number of charge cycles and may eventually need to be replaced.

Fig. e

Should you need to replace the rechargeable battery, remove the small rubber plug from the ACTIVATOR's back cover (Fig. e), unscrew the back cover and remove it. Replace the existing battery pack with a new battery pack (Fig. f) and replace the battery cover. You can purchase a new battery pack by contacting the Careline.

Fig. f

FREQUENTLY ASKED QUESTIONS

The stimulation is uncomfortable. How can I improve this?

- Make sure the GelPads are positioned correctly and that they are pressed firmly against your skin (see pages 7 & 8). Switch off your ACTIVATOR and reposition the garments if necessary.
- You can also smear a few drops of water on the black surface of the GelPads. This can improve the comfort of the muscle stimulation, but be careful not to get water on the ACTIVATOR. Ensure the ACTIVATOR is switched OFF before you do this!
- Ensure the metal studs are fully covered by the GelPads.

My skin is red after the exercise. Is this a problem?

- Some redness of the skin after a toning session is normal. It is partly due to an increase in the blood flow under the skin and should fade after a while. You may also experience some reddening of the skin due to the pressure of the garments. This is the same as the pressure marks you can get from tight clothing. You should not be concerned about this. It should fade soon after you remove them.
- If the redness is excessive, you may have the intensity too high. This may increase the reddening in sensitive skin. Try using a lower intensity for a few days. If the problem persists, you should stop using the product.

Will the product cause muscle soreness?

- As with all exercise some muscle soreness can occur after using ARM LIFT. This is normal and should go away after a while. Use a low intensity for a few sessions if you are experiencing muscle soreness.

How do I know when to replace the GelPads?

- With time the GelPads pick up skin debris and may need to be replaced as this makes the workout less effective and less comfortable.
- You may notice the signal is weakening even if the batteries are OK. This usually indicates that the GelPads are wearing and will soon need replacing. You can check this by first adjusting the garment to ensure correct positioning of the GelPads. If the problem persists, try using the new GelPads. If the signal is still weak or uncomfortable, you should order a new battery pack.
- New GelPads can be purchased online at www.bmrbeautyseries.com or from the Careline at 1-877-518-2323.

Battery performance has degraded significantly.

- qAfter an extended period of time you may notice a degradation in the performance of the ACTIVATOR's battery. At this point you should purchase a new rechargeable battery pack. New battery packs can be purchased through the Careline, or by visiting www.bmrbeautyseries.com.

QUICK GUIDE TO SOLVING PROBLEMS

Problem	Possible Cause	Solution
No signal/ display is blank	The battery is low	Recharge the battery
The display is on but there is no signal	The battery is low	Recharge the battery
The battery symbol is flashing	The battery is low	Recharge the battery
The  symbol has appeared on the display*	The ACTIVATOR and garments are not properly connected	Ensure the garments and ACTIVATOR are properly connected
	Poor GelPad contact with the skin	Press the GelPads firmly against the skin
	The covers are still on the GelPads	Remove the covers from the GelPads
	The GelPads are worn	Visit the website or your local retailer for new GelPads
	The GelPads are not covering the studs	Make sure the GelPads fully cover the metal studs
Unpleasant feeling beneath the GelPads	The GelPads are worn	Visit the website or your local retailer for new GelPads
	GelPads are not covering the metal studs	Reposition the GelPads so that the metal studs are covered
	Too many consecutive sessions	You should just do one session per day one any single body area
Contractions are weak even with a high intensity	The GelPads are worn	Visit the website or your local retailer for new GelPads
	The battery is low	Recharge the battery
	Poor GelPad positioning	See page 7 for correct GelPad positioning
	GelPads are not covering the studs	Reposition the GelPads so that the metal studs are covered
“Err” appears on the display	Error message	Switch the ACTIVATOR off and on again. If the error persists call the Careline

* The  symbol always appears in conjunction with the left/ right intensity indicators. If the left indicator appears, there is a problem with the left-arm garment, the right indicator indicates a problem with the right-arm garment. If both indicators appear there is a problem with both garments.

Caring for your ACTIVATOR

YOUR ACTIVATOR should not be allowed to get wet or be left in excessive sunlight. It may be cleaned regularly using a soft cloth, lightly dampened in soapy water. Do not allow the interior of the ACTIVATOR to become wet. Do not use detergents, alcohol, spray aerosols or strong solvents on the ACTIVATOR.

Access to the interior of the ACTIVATOR is not required for maintenance purposes.

If your ACTIVATOR is damaged, you should not use it, but should return it to BMR or your local retailer for replacement or repair. Repairs, service and modifications may not be carried out by anyone other than qualified service personnel authorized by BMR.

IMPORTANT: Under no circumstance should anything other than the correct type of batteries - (rechargeable batteries 3.6V NiMH) be used with your ACTIVATOR. These can be purchased from the Careline.

Disposing of the GelPads and batteries

Used GelPads and batteries must never be disposed of in a fire but in accordance with national or state laws governing the disposal of such items.

NOTE: The signal gradually increases to a peak intensity level at the start of the contraction phase (i.e. ramp up) and gradually decreases to nothing at the start of the relaxation phase (i.e. ramp down). When the  symbol is flashing, the stimulus is significantly reduced.

Beauty Series Accessories:

You can purchase all accessories at www.bmrbeautyseries.com. Under no circumstances should anything other than BMR LIFT accessories be used with your ACTIVATOR FOR LIFT (Type 390) or charger. Any others may not be compatible with the product and could degrade the minimum safety levels.

BMR TUMMY LIFT Garment (E:70)

- BMR TUMMY LIFT GelPads:**
- 1 large GelPad (Type 706 or 709)
 - 2 small GelPads (Type 707 or 710)

BMR BOTTOM LIFT Garment (E:30)

- BMR BOTTOM LIFT GelPads:**
- 4 large GelPads (Type 716)

BMR ARM LIFT Garments (E:60)

- BMR ARM LIFT GelPads:**
- 4 small GelPads (Type 715)

Belt Extender

Battery Pack (3.6V, NiMH)

Battery Charger (EU - 2504-0303)

Battery Charger (US/Japan - 2504-0302)

Charger complies to IEC 60950 and UL 1950, Input: nominal voltage 100-240V, frequency 50-60Hz, current 90mA. Output: 6V / 550mA. Do not use any other chargers or any other power supplies.

Product Type: 394

Intended use: Muscle stimulator

Classification: Internally powered equipment, Class II charger, Type BF applied parts.

Waveform: Symmetrical bi-phasic square waveform when measured into a resistive load.

Environmental Specifications:

Operating: Temperature Range: 0 - 35°C (32 - 95°F)
Humidity: 20 - 65 % RH

Storage: Temperature Range: 0 - 35°C (32 - 95°F)
Humidity: 20 - 85 % RH

Garment Materials:

Outer material: 80% Nylon, 20% polyester; Inner Material: 100% Polyester laminated to 100% polyurethane; Binding: 100% Nylon; Hook & Loop: 100% Nylon; Foam: EVA; Elasticated Straps: 100% Nylon.

Description of your ACTIVATOR's symbols:

There are a number of technical markings on your ACTIVATOR. These can be explained as follows:

The ACTIVATOR and garments are manufactured for Bio-Medical Research Ltd, Parkmore Business Park West, Galway, Ireland.

The ACTIVATOR requires 1 x 3.6 volt (NiMH) DC battery pack DC is indicated by the symbol: 

The Output Frequency indicates the number of pulses per second transmitted by the ACTIVATOR. This is measured in hertz, which is indicated by "Hz".

Output (RMSA) means the maximum output root mean square current for each channel.

Output (RMSV) means the maximum output root mean square voltage for each channel.

 This symbol means "Attention, consult the accompanying documents".

 This symbol means type BF applied parts.

 This symbol on your ACTIVATOR is to indicate conformity to the requirements of the Medical Device Directive (93/ 42/EEC). 0366 is the number of the notified body (VDE).

 Conforms to UL STD 60601-1.
Certified to CAN/CSA.
STD C22.2 NO 601.1

SN stands for 'serial number'. Inside the box, on the back of the ACTIVATOR is the serial number specific to this ACTIVATOR. The letter preceding the number indicates the year of manufacture, where "P" denotes 2010, "Q" denotes 2011 etc.

The garments' batch number is represented on their packaging, by the number corresponding with the LOT symbol.

Rated Outputs – Voltage/Currents: E:60

Parameter	500Ω	1KΩ	1K5Ω
Output RMSV	6.4V	12.7V	12.7V
Output RMSA	12.7mA	12.7mA	8.5mA
Output Frequency	50Hz	50Hz	50-Hz
DC Component: Approx.	0 C	0 C	0 C
Positive Pulse Width:	175-200µs	175-200µs	175-200µs
Negative Pulse Width:	175-200µs	175-200µs	175-200µs
Interphase Interval:	100 µs	100 µs	100 µs

Attention:

The ACTIVATOR can deliver currents in excess of 2mA/cm² if used with incorrect GelPads.

 At the end of the product lifecycle, do not throw this product into the normal household garbage, but bring it to a collection point for the recycling of electronic equipment.

Some product materials can be re-used if you bring them to a recycling point. By re-using some parts or raw materials from used products you make an important contribution to the protection of the environment. Please contact your local authorities if you need more information about collection points in your area.

Waste Electrical and Electronic Equipment can have potentially harmful effects on the environment. Incorrect disposal can cause harmful toxins to build up in the air, water and soil and can be harmful to human health.

Dear Customer,

Thank you for purchasing this BMR product ("Product"). To be covered by this Limited Warranty, the Product must have been purchased directly from SDI or its authorized retailer in the United States. This Limited Warranty extends only to the original purchaser and cannot be transferred to anyone. Please read this carefully before using the Product.

Our Limited Warranty to you

SDI WARRANTS TO YOU, THE ORIGINAL PURCHASER, THAT UNDER NORMAL USE AND SERVICE, THE PRODUCT WILL BE FREE FROM DEFECTS IN MATERIAL AND WORKMANSHIP FOR A PERIOD OF TWO (2) YEARS FROM THE DATE OF PURCHASE ("WARRANTY PERIOD"). IF YOU DISCOVER A DEFECT COVERED BY THIS LIMITED WARRANTY DURING THE WARRANTY PERIOD, AND RETURN THE PRODUCT TO US, WE WILL REPAIR OR REPLACE THE PRODUCT, AT OUR OPTION, USING NEW OR REMANUFACTURED COMPONENTS. THERE IS NO WARRANTY OR CONDITION OF ANY KIND AS TO ANY DEFECTS DISCOVERED AFTER THE WARRANTY PERIOD EXPIRES.

Certain Product failures not covered

This Limited Warranty covers defects in manufacture that arise from correct use of the Product per the product manual and instructions. It does not cover damage caused by abuse, misuse, accidents, modification or repair by any person other than SDI, moisture, extreme heat or cold, corrosive environments, high voltage surges from external sources or an unsuitable power supply. This Limited Warranty does not cover normal wear and tear on consumable items (e.g. GelPads, garments, etc.). This Limited Warranty does not apply to any battery provided with the Product, or if the Product has an altered or defaced serial number.

How to obtain service under this Limited Warranty

To obtain service under this Limited Warranty, please contact the Customer Careline at 1-877-518-2323 to obtain return instructions. You must submit proof of purchase in the form of a bona fide, dated bill of sale, receipt, or invoice (or copy thereof) evidencing that your request for service is made within the Warranty Period.

Except as may otherwise be required by legislation in your State, SDI will not be responsible for the cost of returning the Product or for any damage to or loss of the Product while it is in transit to us. We recommend that you use a carrier that will track your shipment and ensure the Product's safe delivery. SDI cannot accept any returns sent C.O.D.

THIS LIMITED WARRANTY IS THE ONLY WARRANTY FOR THE PRODUCT, AND THERE ARE NO OTHER EXPRESS WARRANTIES, ORAL OR WRITTEN, PROVIDED BY SDI. ANY IMPLIED WARRANTIES APPLICABLE TO THE PRODUCT, INCLUDING ANY WARRANTIES OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE, ARE LIMITED TO THE DURATION OF THE WARRANTY PERIOD SET FORTH ABOVE.

Limitation of Remedies

YOUR EXCLUSIVE REMEDY, IN LIEU OF ALL INDIRECT, INCIDENTAL, CONSEQUENTIAL, EXEMPLARY, SPECIAL OR OTHER DAMAGES OF ANY KIND, INCLUDING FOR NEGLIGENCE, IS LIMITED TO THE REPAIR OR REPLACEMENT OF THE PRODUCT, IF DEFECTIVE, UNDER THE TERMS AND CONDITIONS OF THIS LIMITED WARRANTY.

SDI will bear no other expenses.

SDI reserves the right to refuse to repair or replace the Product if it is returned in an unsanitary condition. SDI reserves the right to substitute a product of equal or better quality if an identical product is not available at the time of replacement. This Limited Warranty is non-transferable.

Applicability of state law

Some states do not allow limitations on the duration of implied warranties, or the exclusion or limitation of incidental or consequential damages, and so the related limitations and exclusions stated above may not apply to you.

Governing law

This Limited Warranty is valid in the United States, Canada and Mexico only and is not applicable in other countries. If you acquired the Product in the United States, the laws of the State of New York will apply to this Limited Warranty.

THIS LIMITED WARRANTY GIVES YOU SPECIFIC LEGAL RIGHTS AND YOU MAY ALSO HAVE OTHER RIGHTS, WHICH VARY FROM STATE TO STATE.

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