SAD

INSTRUCTION MANUAL
MANUEL D'UTILISATION
BENUTZERHINWEISE
MANUALE DI ISTRUZIONI
GEBRUIKSAANWIJZING
BETJENINGSVEJLEDNING

QUALITY AND RELIABILITY

Here at S.A.D. Lightbox Company we consider quality and reliability our first priority. All components used in the manufacture of our lightboxes have been selected for their proven reliability and we believe that this is reflected in our finished product.

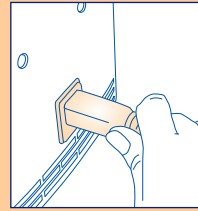
Figure 1

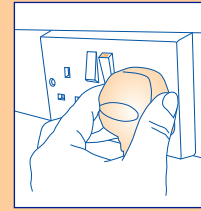
Figure 1a

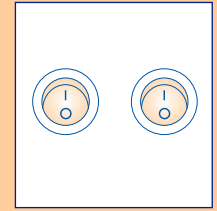
Figure 2

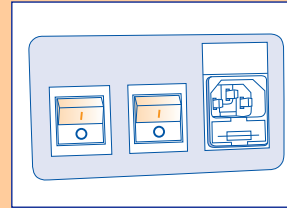
Figure 2a

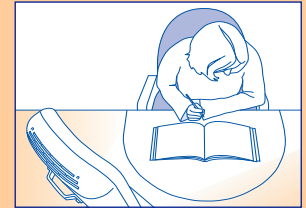
Figure 3

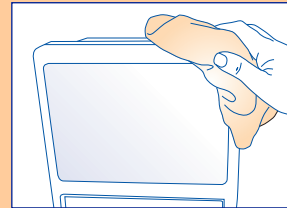
Figure 4

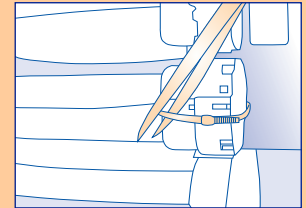
Figure 5

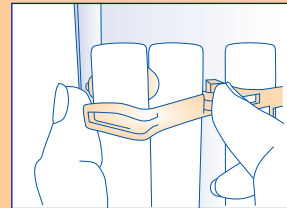
Figure 6

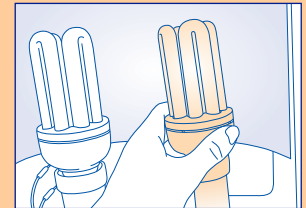
Figure 7

ENGLISH INSTRUCTIONS

S.A.D. is characterised by symptoms such as:

- ▶ Depression, feelings of gloom and despondency for no apparent reason
- ▶ Lethargy, lack of energy
- ▶ Anxiety, inability to cope
- ▶ Social problems, irritability, not wanting to see people
- ▶ Sleep problems, finding it hard to stay awake during the day, and having disturbed nights
- ▶ Loss of libido, not interested in sex or physical contact
- ▶ Craving for carbohydrates and sweet foods, increased appetite leading to weight gain

It generally occurs during the darker winter months of September through to March or April. Health experts agree that the winter blues respond - often dramatically - to bright light therapy. Some people do use their lightboxes all year, but for the majority it is an autumn/winter treatment. The important thing is, do not wait for the Winter Blues to take hold, start using your lightbox before the clocks change, say in late August or September and get into the rhythm of light treatment.

PREPARING THE LIGHTBOX FOR USE

- 1 Place the lightbox on a firm surface, ideally white or light coloured.
- 2 Plug the cable into the back of the lightbox and into the wall socket.
Figure 1/1a.
- 3 Switch the lightbox on. The switches will illuminate. *Figure 2/2a.*

DIAMOND LIGHTBOX INSTRUCTION FOR FITTING BASE, FLOOR STAND AND BRACKET

BASE

Parts included:

- Base x 1
- Pressure plate x1
- Thumb screws x 2
- Thumb screw retainers x 2
(these are white, and only secure the parts together for transit)

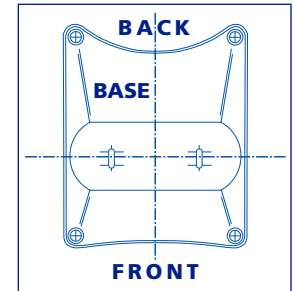
Diagram 1

- 1 Remove the lightbox and base from carton.
- 2 Place lightbox on its back on a table or firm surface.
- 3 Place one hand under assembly to hold the base and the pressure plate together, then remove the small white thumb screw retainers before step 4.
- 4 Offer up the base to the lightbox which fits neatly into the depression in the base, and the holes that the screws need to be inserted into will all line up.
- 5 The thumb screws just need to be tightened up as much as possible to secure the lightbox to the base.

*N.B. If the lightbox faces the front of the base as seen in **Diagram 1**, then the angle of tilt forward will be 6° forward and 12° backwards. Should you wish to, you can secure the lightbox to the stand the other way round to increase the forward tilt to 12°, and the backward tilt to 6° depending on which angle is most convenient for you.*

OPTIONAL EXTRAS FOR DIAMOND RANGE

FLOOR STAND

Parts Included:

- Bracket x 1
- Floor stand x 1
- Thumb screws x 2
- Thumb nut x 1
- Thumb screw retainers x 2
(these are white, and only secure the parts together for transit)


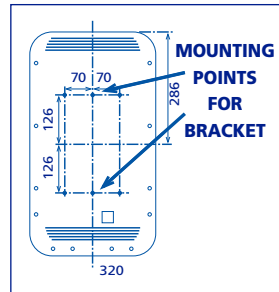
Diagram 2

Tools required: A flat blade screwdriver

- 1 Remove handle and base from lightbox if already fitted.
- 2 Then remove the small white retainers from the bracket.
- 3 Place lightbox face down on a table or firm surface.
- 4 Offer the bracket up to the two central mounting points as shown in **Diagram 2**, and tighten thumb screws. *N.B. The larger of the holes in the bracket goes towards the base of the lightbox.*
- 5 Set up the stand by extending the central pole to the maximum length and lock the bottom knob to hold secure. (The pole is telescopic. To release piece inside you need to undo the long 'collar' and extend the rest of the pole. Tighten the collar again.)
- 6 Remove the thumb nut.
- 7 Lift lightbox up and slot holes over the pole of the floor stand and push down until you are able to replace thumb nut, and tighten.

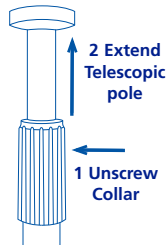WALL BRACKET

Parts Included:

- Bracket x 1
- Thumb screws x 2
- Thumb nut x 2
- Thumb screw retainers x 2 (these are white and only secure the parts together for transit)

Items required: A screwdriver

- 4 Screws (we recommend size 8 x 1 1/2" (38mm) minimum)
- 4 Rawl plugs for the type of wall in your home

- 1 Remove handle and base from lightbox if already fitted.
- 2 Remove the small white thumb screw retainers from bracket.
- 3 Place lightbox face down on a table or firm surface.
- 4 Offer the bracket up to the two central mounting points as shown in **Diagram 2**, with the rubber grommet at the bottom of the lightbox, and tighten thumb screws.
- 5 The wall bracket then needs to be attached to the wall.
- 6 Remove the two centre thumb nuts and lift the lightbox up and offer it up to the wall bracket. The wall bracket can be attached to the lightbox bracket at either the top, middle or bottom fixing points as required.
- 7 Tighten the thumb nuts, adjust position of lightbox and then tighten the T bar to lock the universal ball joint.

TAKING YOUR TREATMENT

- 1 Have the lightbox about an arms distance away and slightly to one side. **Figure 3.**
- 2 The light needs to go into the eye, so do not wear goggles or tinted glasses. Plain glasses or contact lenses are fine. Do not doze or sleep - you must be awake to take the treatment.
- 3 Do not look at the light all the time, just glance at it every now and then. You can carry on with your normal routine, reading and writing; have breakfast beside it or use when you are using your computer.
- 4 **It is better to take your light treatment as early as possible during the day, and you should not use it within 4 hours of bedtime.**
- 5 Repeat daily during the autumn and winter months.

LENGTH OF TREATMENT

Lightbox	Replacement Tubes <i>(See below for instructions on how to change tubes)</i>	Daily Treatment Time
Travel Lite	ND23W	Approx 1¼ hrs
Diamond 2, Hi Lux	ND36W	Approx 1 hr
Diamond 3		
Super Lux, Super Lite 3	ND36W	Approx 30 mins
Diamond 4, Silver Lite	ND36W	Approx 20 mins
Ultima	ND40W	Approx 20 mins

These are average times and if your symptoms are severe then times can be increased. If you are working you can use your lightbox for longer than recommended but not within 4 hours of bedtime.

SAFETY

DO NOT

- › Use near water or in damp conditions
- › Use in the presence of flammable gases
- › Cover the lightbox with anything that could obstruct ventilation slots

MAINTAINING YOUR LIGHTBOX

Always switch the lightbox off at the switch and remove plug before maintenance.

CLEANING

Clean occasionally with a soft duster or wipe with a damp cloth **figure 4.**

CHANGING TUBES

We recommend that tubes are changed every 3 years to ensure a high light output.

TO CHANGE LAMPS IN THE WOODEN AND PLASTIC RANGE

To change lamps (this should only be undertaken by a competent person or electrician.):

- 1 Disconnect lightbox from the mains supply and remove screws to take box apart (on the side of the plastic box, on the back of the wooden box*).

- 2 Cut the retaining cable tie at the bottom of the light and discard (replacement of this item is not necessary as it is used to retain the lamp during transit). **Figure 5.**
- 3 Remove the plastic lamp clip at top of lamp (this is a two part moulding and can be removed by placing a thumb under one side and flicking it free, it can sometimes be quite stiff), and slide lamp out. **Figure 6.**
- 4 Insert new lamp and replace lamp clip.
- 5 Reassemble the box, ensuring that the reflector is correctly positioned and replace screws. (In the wooden boxes there are 6 small plastic clips around the screen where the reflector locates. If it does not, you can carefully insert a fine screw driver through slots in the front of the diffuser and 'flick' it into place.)

** If the back of the wooden box will not come out, insert some longer screws into the screw holes. Leave these proud and they will give you leverage to remove the back which can sometimes be a little stiff.*

▶ TO CHANGE LAMPS IN THE DIAMOND RANGE

As above, but the screws are located on the back of the box. Only the screws holding the diffuser need to be removed. (Do not remove the bottom 2 screws as this removes the ballast cover.)

▶ TO CHANGE LAMPS IN THE TRAVEL LITE

- 1 Disconnect lightbox from the mains supply and remove the screws around the edge of the box, and remove the front.

- 2 Unscrew the old lamp by holding the base not the tubes, these can be quite tight. **Figure 7.**
- 3 Insert and screw in the new lamp and tighten as much as possible. It needs to be screwed down as far as possible to make a good connection.
- 4 Replace front.

N.B. The two longer screws go at the top of the light box.

TECHNICAL SPECIFICATION

220-240Vac 50Hz
Operating temperature -10-50°C
Safety EN60601
Class IIa Medical Device

GUARANTEE

All lightboxes are guaranteed for 3 years (tubes 1 year) against manufacturing defect from the date of purchase. For warranty returns see back of instructin manual.

MANUFACTURER:

FULL SPECTRUM LIGHTING T/A S.A.D. LIGHTBOX COMPANY

Unit 2, Aston Hill, Lewknor, Watlington, Oxfordshire
OX49 5SG • UK
TEL 0845 095 6477
FAX 0845 095 6499
EMAIL info@sad.uk.com

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