

Peak Performer 1 TCD (940- Series) User Manual



Peak Laboratories, LLC

www.peaklaboratories.com

650-691-1267

Version 5.0

1/30/14

Forward

This manual is a reference guide for the operation of the Peak Performer 1 series of gas analyzers. The Peak Performer 1 should only be operated by trained personnel familiar with the safe operating procedures of the analyzer.

The following cautions should always be considered:



High voltages are required for proper detector operation and are present when the unit is energized.



High operating temperatures are required for proper detector and chromatographic operation. Hot surfaces may be present when the unit is energized.



Illustrations and photographs contained herein are for example only. Actual configurations may vary depending on specific application and installed options.



Improper installation, operation or service of the analyzer can cause permanent damage to the instrument. The analyzer is designed to be operated with all covers installed, connected to a properly grounded 3 conductor AC line source.



In event of fire, Class A, B, or C fire extinguishers can be used.



The Peak Performer 1 must have standard electrical power and applicable gas supply pressures as noted on the test documentation for proper operation. Permanent damage and voiding of the warranty may result from improper use.



Operators should not attempt to repair the instrument except under directed to do so by factory trained service technicians. Permanent damage and voiding of the warranty may result of improper operation.

Glossary and Terms

TCD	Thermal Conductivity Detector
Molecular Sieve	Zeolite column packing material
Unibeads	Silica column packing material
HSD	HayeSep D column packing material
Unibeads	Silica column packing material
Gas Purity	Amount of undesired elements in a gas supply
Span Gas	Gas supply with known, certified amounts of specific compounds
Response Factor	Proportionality factor between area count units and concentration units
Name	Peak identification tag
PkCen	Expected peak retention time (in seconds)
LW	Typical time span measured from the start of the peak's baseline rise to the peak apex " <i>PkCen</i> " (in seconds)
RW	Typical time span measured from the peak's apex " <i>PkCen</i> " to the end of the peak's baseline decline (in seconds)
PkWin	Total tolerance window (in seconds) for assignment of a " <i>Name</i> " to a quantified chromatographic peak. The tolerance window is centered upon the " <i>PkCen</i> " value
PkHgt	Variable for establishing the cross-over point in peak detection modes. Chromatographic peaks higher than this parameter will be quantified using the " <i>Variable</i> " mode heights, peaks with height less than or equal to " <i>PkHgt</i> " will be quantified using the " <i>ForceB</i> " mode.
Flt	Convolution filter value expressing the overall peak shape. Flt = 2 is recommended for sharp, narrow peaks such as H ₂ , and Flt = 8 is recommended for broad peaks.
VICI	Valco Instrument Company, Inc. (www.vici.com)

Table of Contents

1.0	Introduction	7
2.0	Installation Considerations	8
2.1	Unpacking Instructions	8
2.2	General Considerations	8
2.3	Electrical Connections	8
2.4	Rack Mount Specifications	9
2.5	Gas Supplies and Connections	11
2.5.1	Carrier Gas Supply	12
2.5.1.1	Carrier Gas Purity	12
2.5.1.2	Carrier Gas Connection	13
2.5.2	Actuator Gas Supply	13
2.5.3	Sample Supply	14
2.5.3.1	External Calibration Gas Supply	14
2.5.3.2	Second Sample Gas Supply	14
2.5.4	Calibration Gas Connection for Analyzer with Internal Blender	15
2.5.4.1	Span Gas Supply (For Internal Span Gas Blender)	15
3.0	Standard Start-up Procedure	16
3.1	Quick Start Sequence	16
3.2	Column Conditioning	18
3.3	Detector Conditioning	19
3.4	Sample Analysis	19
3.5	Calibration	19
3.5.1	Establishing New Instrument Response Factors	19
3.6	Shutdown and Transport	20
4.0	Peak Performer 1 Operator Interface	21
4.1	Set-up Screen	21
4.1.1	Set-up Screen Upper Section	22
4.1.1.1	Set-up Screen Upper Section Button Functions	22
4.1.2	Set-up Screen Middle Section (Power-up)	23
4.1.2.1	Set-up Screen Middle Section Button Functions (Power-up)	23
4.1.3	Set-up Screen Middle Section (Normal)	24
4.1.3.1	Set-up Screen Middle Section Button Functions (Normal)	24
4.1.4	Set-up Screen Lower Section	25
4.1.4.1	Error Messages on the Set-up Screen Lower Section	25
4.1.4.2	Clearing Error Messages on the Set-up Screen	26
4.1.4.2.1	Set-up Screen Lower Section Button Functions	26
4.2	Analog Interface (Trend and Recorder Output) Screen	28
4.2.1	Analog Interface Screen Button Function	28
4.2.2	Analog Interface Connections	29
4.2.3	Analog Interface Wiring Breakdown	30
4.3	Calibration Database Screen	31
4.3.1	Calibration Screen Button Functions	31
4.3.2	Calibration	31
4.3.2.1	Establishing New instrument Response Factors	32
4.4	Analysis Database Screen	33
4.4.1	Analysis Screen Parameters	33
4.4.2	Analysis Screen Button Function	34

4.5	Event Program Editor Screen	35
4.5.1	Event Program Commands	35
4.5.2	Event Editor Screen Button Function	36
4.6	Stream Selector Sequence Editor Screen	36
4.6.1	Stream Selector Commands	36
4.6.2	Stream Selector Program Editor Screen Button Function	37
4.7	Run / Chromatogram Screen	37
4.7.1	Run Screen Button Functions	38
4.7.2	Run Screen / Continuous Monitoring Analysis	39
4.7.3	Manual Screen (Sub-Screen of Run Screen)	39
4.7.3.1	Manual Screen Button Functions	40
4.7.4	Display Screen (Sub-screen of Run Screen)	40
4.7.4.1	Display Screen Button Functions	40
5.0	Viewer Software	41
5.1	System Requirements	41
5.2	Installing the Viewer Software	41
5.3	Starting the Viewer Software	41
5.4	COM1 / COM2 Port and Cable Wiring	41
5.5	Viewer Software Window	42
5.5.1	File Menu on the Viewer Taskbar	43
5.5.1.1	Communication Setting Window	44
5.5.1.2	Viewer Line Details Window	46
5.5.2	Start Menu on the Viewer Taskbar	46
5.5.3	Stop Menu on the Viewer Taskbar	47
5.5.4	Chrom View Window and Chromatogram Viewer Window	47
5.5.4.1	Chromatogram Files	48
5.5.4.2	File Menu on the Chrom Viewer Taskbar	48
5.5.4.3	Viewer Menu on the Chrom Viewer Taskbar	49
6.0	Communication Port Protocols and Data Formats	50
6.1	COM1 / COM2 Port and Cable Wiring	50
6.2	Port 1 Protocol	50
6.2.1	COM 1 Port Settings	50
6.2.2	Peak View Output (Columns format in .CSV file)	50
6.2.3	Detector Output	51
6.2.3.1	External Start Command String	51
6.2.3.2	Detector Signal Data String	51
6.3	Port 2 Protocol	52
6.3.1	COM 2 Settings	52
6.3.2	COM 2 Port PLC Output Format	53
6.4	Modbus Protocol and Register Map for Peak Devices	53
6.4.1	Mode of Transmission	53
6.4.2	Modbus Packet Structure	54
6.4.3	Packet Communications	55
6.4.4	Modbus Data Format	56
6.4.5	Broadcast Packets	58
6.5	Trend / Analog Outputs	60
6.6	Setting up a RS485 Network	60
6.6.1	RS232 to RS485 Converters	60

7.0	Chromatographic Principles	61
7.1	Sample Injection	61
7.2	Carrier Gas Purification	61
7.3	Component Separation	62
7.4	Flow Venting and Balancing	64
7.5	Effects of Temperature and Flow on Analytical Columns	67
7.6	Peak Identification and Quantification	69
7.6.1	Peak Quantification: Fixed Mode	73
7.6.2	Peak Quantification: Forced Baseline Mode	73
7.6.3	Peak Quantification: Variable Mode	74
8.0	Thermal Conductivity Detector	78
8.1	Principles of Operation of Thermal Conductivity Detection	79
9.0	Analyzer and Maintenance and Service Procedures	82
9.1	Peak Performer 1 System Block Diagram	83
9.2	Monitor Regulator Pressures	84
9.3	Verifying Valve Integrity (Carrier Blank, Zero Gas)	85
9.4	Adjusting the Display Contrast	86
9.5	Altering the Analytical Range of the Instrument	87
9.6	Testing the Heater Resistance	88
9.7	Main DC Power Supply Replacement	89
9.8	Adjusting the 5 Volt DC Supply	91
9.9	COM 1 Connections on Multiple Analyzers to Signal Computer	92



1.0 Introduction

The Peak Laboratories Peak Performer 1 (PP1) TCD gas analyzer is a trace level gas analysis system capable of detecting impurities in a variety of sample matrices.

Due to the specific nature of the detection method, analysis times are rapid and chromatographic complexity is minimized. The result is a simple and reliable system capable of following rapid trends in component concentration in laboratory, field, and process gas applications.

The Peak Performer 1 has a fully integrated, stand-alone microprocessor operating system. Key functions of the operating systems include detector and temperature control, data collection and analysis, and operator / data interface

Sampling hardware	4, 6, and 10 port VICI valves
Column oven	Covered
Type	Isothermal, Mandrel Heating
Capacity	Single or Multiple Columns
Temperature range (Ambient)	15°C → 30°C
Temperature stability	± 0.1°C
Power required	160 watts maximum

2.0 Installation Considerations

2.1 Unpacking Instructions:

Remove the Peak Performer 1 from the packing material or traveling case, taking care to not scar the exterior surfaces of the analyzer.

Set the analyzer on a firm, even surface and remove the protective plastic wrap.

Each analyzer ships with adequate supplies for installation. Be certain to remove the installation accessories (tubing, power cord, compression fittings, etc.) from the packing materials.

Remove the sleeve containing the Certificate of Conformity and Installation CD from the analyzer top panel.

All shipping materials are reusable, please recycle these materials appropriately.

2.2 General Considerations

- Analyzer size is 26" L x 17" W x 7" H
- Maximum operating environment for the PP1 analyzer is 30°C
- Minimum operating environment for the PP1 analyzer is 15°C
- Power consumption is 160 watts maximum.
- Fuse size is 2.5 ampere @ 250 VAC, 5 x 20 mm SLO-BLO
- Free air flow for adequate ventilation to the rear of the unit is required.
- Communication between the PP1 and your PC is via 9-pin straight through cable connection to COM 1 and COM 2.
- Analog outputs are available from the DIN screw terminal strips on the rear panel

2.3 Electrical Connections

Verify the correct operating voltage as marked on the rear of the analyzer at the main power switch.

Attach a standard power cable to the rear of the instrument (see Figure 1) and connect to an appropriately grounded outlet.



Figure 1. AC Power Connection Details

The PP1 has been CE certified to be immune to AC frequency and voltage variances of +/- 10% of the nominal AC operating voltage. Refer to the final test report and/or markings on the rear of analyzers to verify the proper operating voltage for each unit.

For best stability, power to the analyzer should originate from an electrical circuit free of large inductive or other current loads. Be sure to mount the analyzer in an area of adequate ventilation and make sure the cooling fan inlet is unobstructed.

There are no battery back-ups or reserve power supplies built into the PP1, therefore any disruption of the AC supply will result in restarting of the analyzer. Peak Labs recommends the installation of a high quality on-line or double conversion type of uninterruptible power supply (UPS).

Be sure to mount the analyzer in an area of adequate ventilation and make sure the cooling fan inlet is unobstructed. The operating environment for the PP1 analyzer should be maintained between 15°C and 30°C.

2.4 Rack Mount Specifications

The PP1 is designed for mounting in a standard 19" instrument rack. The analyzer occupies 4U (6.9") of panel height, 26" in depth.

Rack retaining flanges provided with the Peak rack mount kit attach via the analyzer top cover to provide secure connection to the instrument rack rails.

Peak recommends use of high quality 26" full extension ball bearing slides for instrument mounting, such as Jonathan QD375-26 or equivalent.

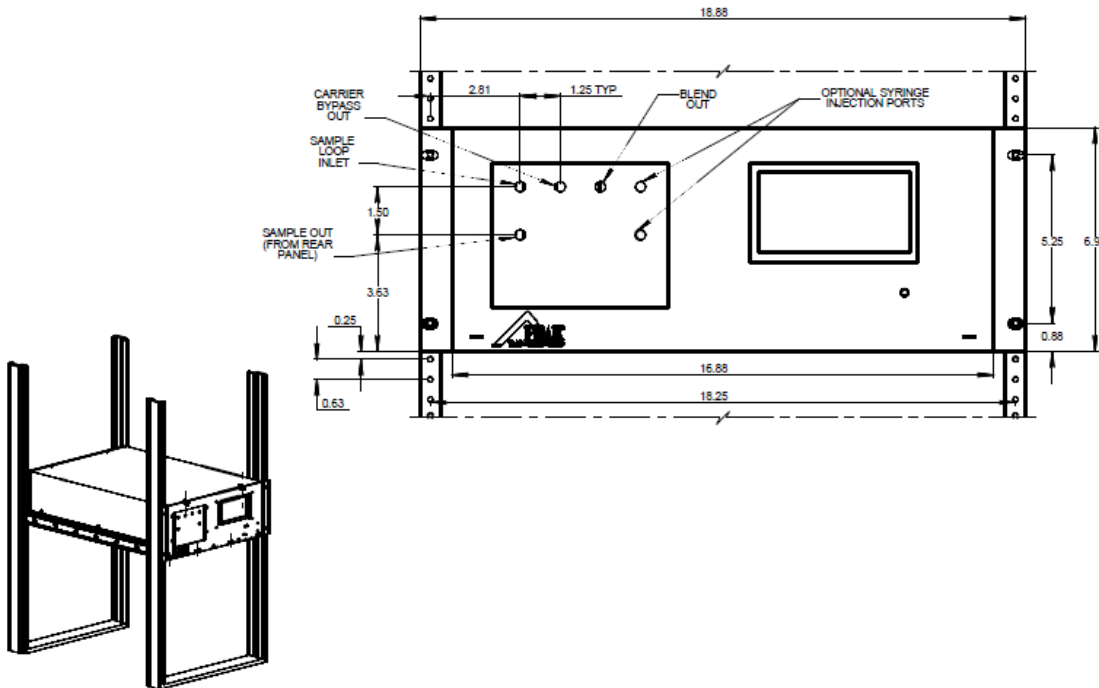


Figure 2. Rack Mount Front View with Dimensions

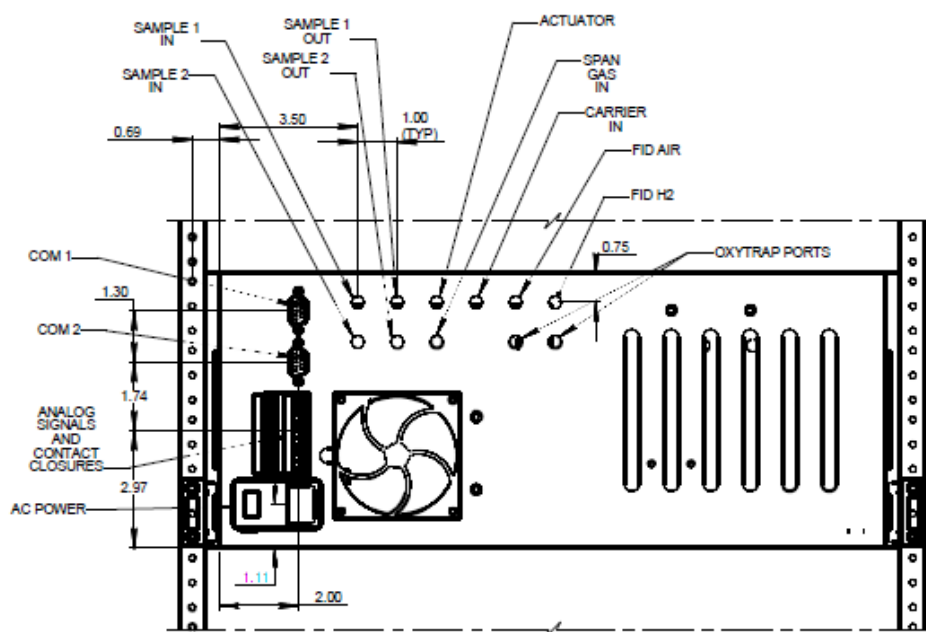


Figure 3. Rack Mount Rear View with Dimensions

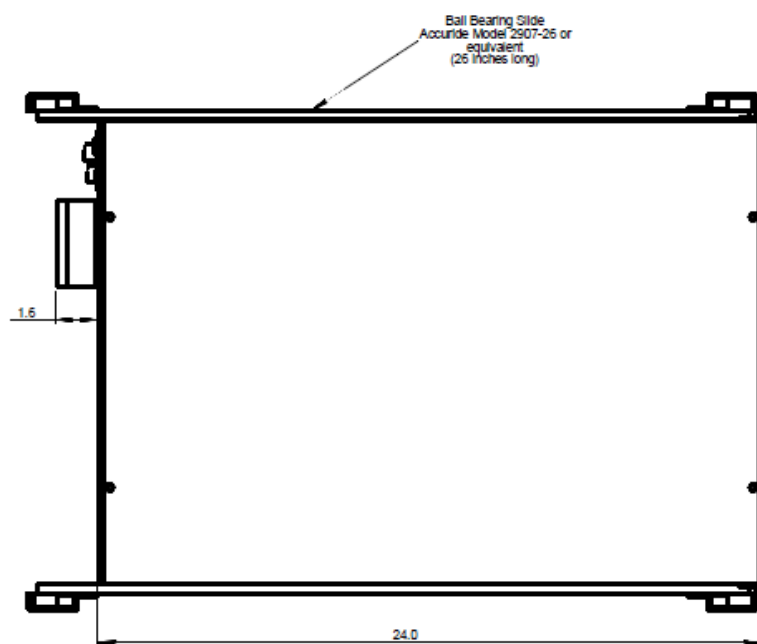


Figure 4. Rack Mount Top View with Dimensions

2.5 Gas Supplies and Connections

All plumbing connections are 1/16" female VICI bulkhead fittings.

The Peak Performer 1 is shipped with tubing and fittings required to connect the analyzer to gas sources.

Replacement materials are available from Peak Laboratories, or may be purchased directly.

- Standard gas fittings are 1/16" VICI compression fittings (VICI P/N ZN1 & ZF1)
- Supply tubing is 1/16" O.D x 0.03" I.D, cleaned and baked T300 stainless steel tubing



Figure 5. Rear View of PP1 (Some Fittings May Not Be Present on Certain Models)

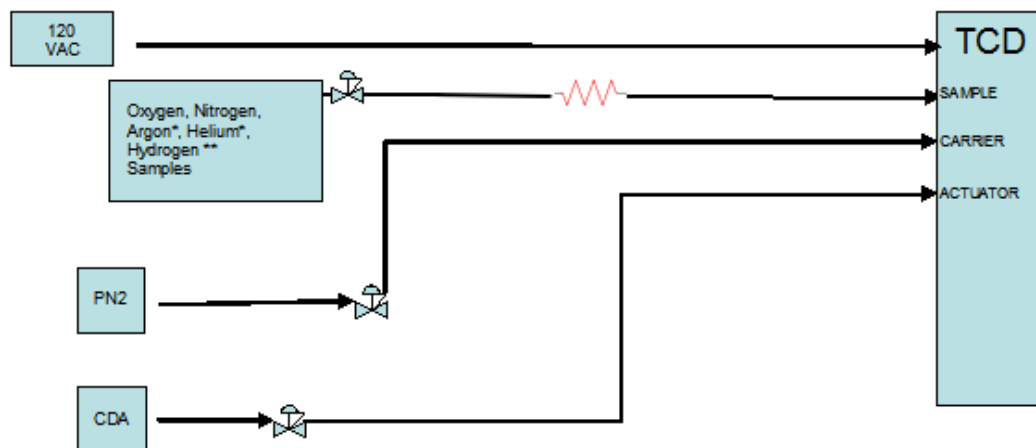


Figure 6. TCD PP1 Gas Supply Schematic

2.5.1 Carrier Gas Supply

The instrument will accept a variety of carrier gases including Ar or He. Please consult a Peak representative if you would like to use other carrier gases.



**** CAUTION**:** Oxygen is never a suitable carrier or actuator supply gas. Oxygen is highly reactive, and many pneumatic components contain greases and oils that combust spontaneously when exposed to oxygen.

Consult the test data shipped with the analyzer (in CD format) for carrier gas details specific for your analyzer.

Commonly, the PP1 will utilize a purified carrier gas supply, and helium is considered the best multi-purpose TCD carrier gas.



PLEASE CONSULT THE TEST DATA ATTACHED TO YOUR ANALYZER IF UNCERTAINTY OF CARRIER GAS SELECTION EXISTS.

2.5.1.1 Carrier Gas Purity

Chromatographic instrument detection limit is directly related to carrier gas purity. Improved carrier gas purity enables improved sensitivity.

Typical Carrier Gas Specifications (Pre-Purification)

Source	High Pressure Cylinder or Liquid Dewar
Gas Type	Ar or He
THC Concentration	< 10 ppm
CO, CO ₂ Concentration	< 10 ppm
O ₂ Concentration	< 10 ppm
H ₂ , H ₂ O Concentration	< 10 ppm

Consequently, the analyzer's minimum detectable quantity (MDQ) would be quite high (> 1000 ppm) unless the carrier gas is purified.

2.5.1.2 Carrier Gas Connection

The PP1 carrier gas connection utilizes a 1/16" female VICI bulkhead labeled "CARRIER IN" on the front or rear panel.

- 1) Connect carrier gas supply gas (Ar or He) to the carrier bulkhead and set the source regulator to the setting shown in the Final Test data.
- 2) Verify detector flow with a flow meter to roughly match the setting in the Final Test data.
- 3) Verify the "BYPASS OUT" is flowing at least 15 sccm.

Carrier flow entering the analyzer is split with a portion of the flow passing through a gas sampling valve and another portion is passing through a restrictor tube terminating at the location of the "BYPASS OUT" port. Under normal operating conditions the "BYPASS OUT" port is plugged.

During column reconditioning, the plug may be removed and the "SAMPLE IN" switched to the bypass flow. In this manner, column reconditioning may be performed without requiring a shutdown of the sample gas flow, avoiding the lengthy re-equilibration period.

2.5.2 Actuator Supply

The actuator supply is normally independent of the carrier gas inlet. Actuator supply is normally attached to the rear panel of the instrument. Attach the air or inert gas actuator supply to "ACTUATOR IN" port on the rear panel. Preset the source to 60 - 80 psig.

Gas Supply	Gas	Purity	Nominal Pressure	Typical Flow
Actuator	Clean, Dry Air or better	< 10 ppm H ₂ O	70 psig	100 sccm

2.5.3 Sample Supply

Sample lines are normally attached to the rear panel of the instrument. Attach sample or span gas lines to the SAMPLE IN port on the rear panel. The SAMPLE OUT port must be unrestricted.

Sample gas must be supplied to the analyzer at low pressure (0.3 to 1.0 psig) and moderate flow (30 to 120 ml/min). The analyzer does not have any sample gas pressure / flow regulation equipment. Higher sample pressures can be used if a fixed restrictor is set and calibrated to the incoming pressure. Please contact the factory for sample handling instructions if these input conditions are not available.

SAMPLE GAS SUPPLY

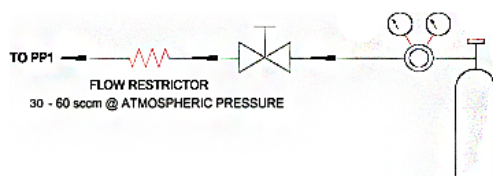


Figure 7. Sample Supply Flow Schematic

2.5.3.1 External Calibration Gas Supply

Certified calibration standards must be connected periodically to the PP1 for verification of analyzer response. When using a mobile calibration gas source such as high pressure cylinder, the plumbing arrangement shown in figure 7 is recommended. The calibration source can be connected to the SAMPLE IN port on the analyzer front panel, or alternatively, make sure the front panel jumper between SAMPLE OUT and SAMPLE IN is installed and connect the calibration source to the rear panel SAMPLE IN port. The SAMPLE OUT port must be unrestricted when using either method.

2.5.3.2 Second Sample Gas Supply

Peak Labs offers an optional second sample inlet stream for most analyzers. Connection to the SAMPLE 1 IN and SAMPLE 2 IN ports should replicate Figure 7. Sample lines are normally attached to the rear panel of the instrument. Attach sample or span gas lines to the SAMPLE IN ports on the rear panel. The SAMPLE OUT ports must be unrestricted.

2.5.4 Calibration Gas Connection for Analyzers Equipped with the Optional Internal Span Gas Blender.

Peak Labs offers an optional internal span gas blender for most analyzers. Connection to the SPAN GAS IN port should replicate Figure 8. Nominal inlet span gas pressure is 50 psig and flow through the internal blender is preset within the analyzer. External restrictors are not required. It is important to measure the pressure applied to the SPAN GAS IN port as closely to the analyzer rear panel as possible.

2.5.4.1 Span Gas Supply (Optional: For internal span gas blender units only)

Peak Labs offers an optional internal span gas blender for most analyzers. Connection to the SPAN GAS IN port should replicate Figure 8. It is important to measure the pressure applied to the SPAN GAS IN port as closely to the analyzer rear panel as possible. Nominal inlet span gas pressure is 50 psig and flow is restricted inside the analyzer. External restrictors are not required for these units

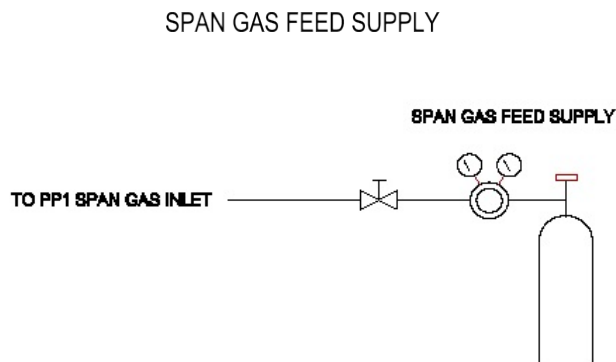


Figure 8. Span Gas Supply Flow Schematic

3.0 Standard Start-up Procedure

Peak Labs recommends that the user read the entire operating manual prior to using the “Quick Start” sequence.



**** CAUTION **** Do not begin detector heating without carrier gas flow. Damage to the detector may result.

3.1 Quick Start Sequence

1. Confirm that all gas connections and supplies are properly made.
2. Remove the top cover of the instrument.
3. Pressurize the carrier gas supply line as specified in the Final Test data.

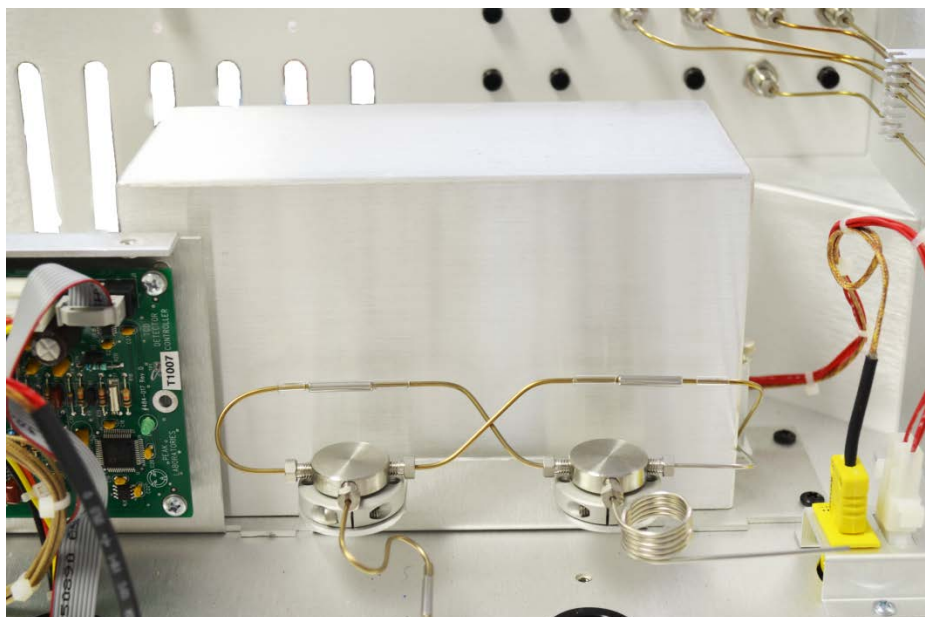


Figure 9. TCD Detector

4. Verify that sample is flowing freely through loop - exit flow from the rear SAMPLE OUT port should be 20 to 120 cc/min.
5. Verify “BYPASS OUT” flow > 20 cc/minute.
6. Verify the analyzer is connected to the correct line voltage supply.
7. Energize the main power switch (located to the right of the power cable, see Figure 1).
8. After energizing, the instrument will display the MAIN screen in the IDLE state.
9. Set Date and Time
10. Press “↓” command three times to position the cursor in the middle 1/3 of the screen.
11. Press the “HEAT” command once to activate heater zones.
12. Allow the 20 -24 hours for the analyzer to thoroughly stabilize the heated zones.
13. After temperature stabilization, energize the filament

14. Verify that the TCD Vout reading is greater than 50000 mV and less than 200000.
15. The TCD Quality number to the far right of the TCD Vout reading should be between 4 and 40
16. If the analyzer has been stored or idle without carrier gas flow for more than a few days, consider column conditioning as described in Section 3.2
17. Press the "PAGE" command 5 times to reach the RUN SCREEN.
18. If automatic sample cycles are desired, press the "SINGLE" command once to change the setting to "CYCLE"
19. To initiate sampling, press the "RUN" command once. The P1 event program will be automatically loaded and run. Cycle samples with "SAMPLE IN" connected to "BYPASS OUT" on the front panel (See Step 17 to start "CYCLE" process)



NOTE: To display concentration data only, press the upper right corner of the RUN screen once – this will hide the chromatogram and display data in large text format. To reveal the chromatogram again, press the upper right corner of the RUN screen again.

20. To halt sampling, press the "IDLE" command once.

3.2 Column Conditioning

Conditions warranting column reconditioning are:

- a) Start-up after storage / shipment
- b) Reduced detector response,
- c) Unstable baseline
- d) Merging of peaks.

Normal column operating temperatures are 45° to 150°C. The maximum PP1 TCD column operating temperatures 250°C.



***** CAUTION ***** Never energize heaters of the instrument unless carrier gas flowing.



***** CAUTION***** Do not exceed 300°C of TCD column oven temperature. Temperatures in excess of 300°C can destroy column material.

- 1) Connect the front jumper tube between SAMPLE IN and BYPASS OUT.
- 2) Use the ↑, ↓, ← and → keys to navigate and alter the SET value for the "Column" heater zone to 180° to 200°C.
- 3) After 8 to 12 hours of stable conditioning temperature, the column temperature may be restored to the original factory setting.
- 4) After normal operating temperature is achieved, the front jumper tube may be reconnected between SAMPLE IN and SAMPLE OUT.
- 5) After a 5-10 minute detector re-equilibration period, samples may be run.

Normal operating temperatures vary with specific application. Consult the Final Test Data for special settings, special reconditioning requirements, etc. for your application.

3.3 Detector Conditioning



**** CAUTION **** The detector filament operating temperature is preset at the factory. The

TCD detector does not normally require conditioning treatment, however it does require 20 to 24 hours at normal operating temperature to achieve optimum stable performance.

Normal detector operating temperature for the TCD is 150°C, although operating temperatures vary with the specific application.

3.4 Sample Analysis

Verify that sample is flowing freely through loop - exit flow from the rear SAMPLE OUT port should be 30 to 120 cc/min.

Press the "ESC" or "PAGE" commands as needed to navigate to the RUN screen

If multiple sample runs are desired, press the "SINGLE" command to toggle to "CYCLE" (repetitive analysis) mode

Press the "RUN" key to initiate analysis – the P1 event program will be loaded automatically and executed. Allow the analyzer to fully execute the event program. Current concentration data will be displayed as an overlay on the RUN screen as peaks are identified and quantified.

3.5 Calibration

The Peak Performer 1 chromatograph determines concentration peak area to quantify the amount of contaminants in the sample gas. The direct instrument response (in units of uV-sec) is reported simultaneously as an overlay on the RUN / CHROMATOGRAM SCREEN and in the AREA column on the CALIBRATION DATABASE SCREEN.

As the concentration peak area is determined, the analyzer applies the pre-programmed response factor as an inverse proportional factor to generate the reported concentration value.

$$\text{PEAK CONCENTRATION} = \text{PEAK AREA RESPONSE} / \text{RESPONSE FACTOR}$$

3.5.1 Establishing New Instrument Response Factors

Connect a known calibration standard gas (span gas) to the analyzer at the SAMPLE IN port on the front panel. (If span gas connection is made at the rear panel SAMPLE IN port, ensure that the front jumper is connected between SAMPLE OUT and SAMPLE IN on the front panel)

Verify that sample is flowing freely through loop - exit flow from the SAMPLE OUT should be 30 to 120 cc/min.

Press the ESC or PAGE keys as needed to navigate to the RUN screen

Press the SINGLE / CYCLE / RERUN key to toggle to SINGLE mode

Press the RUN key to initiate analysis – the P1 event program will be loaded automatically and executed. Allow the analyzer to fully execute the event program. Current concentration data will be displayed as an overlay on the RUN screen as peaks are identified and quantified.

Press the PAGE key as needed to toggle to the CALIBRATION DATA BASE screen

Place the “ * ” at the first compound of interest and press the EDIT key – use the up, down, left and right arrow keys to enter the span gas concentration in the CONC field. Be certain to use the right arrow key to fully exit the CONC field or all updates will be lost.

With the “ * ” before the compound name, press the UPDATE key. The analyzer will recalculate the appropriate response factor for the compound based upon the span gas concentration and analyzer area response using the following formula:

$$\text{RESPONSE FACTOR} = \text{AREA RESPONSE} / \text{SPAN GAS CONCENTRATION}$$

NOTE: The LED on the analyzer's front panel will change from green to red to remind the user that the current calibration information is stored in volatile (temporary) memory only, and will need to be saved to permanent memory at a later time.

Repeat the sequence for each of the compounds shown on the CALIBRATION DATA BASE screen.

Repeat analysis of the span gas sample and verify that the analyzer reports data values within normal calibration tolerances (i.e. + / - 5%).

3.6 Shutdown and Transport

1. Allow the temperature zones to cool below 50°C (approximately 2 hours).
2. Allow carrier to flow during cool down time.
3. Shut off carrier after cool down time and plug the CARRIER IN port on the rear of the analyzer.
4. Power off analyzer.
5. Cap the SAMPLE IN, SAMPLE OUT and ACTUATOR IN ports.
6. Plug the BYPASS OUT port on the front panel.
7. Ensure that the front jumper tube is connected between SAMPLE IN and SAMPLE OUT.
8. Package the analyzer in a hard cover carrying case or comparable shipping container. If cardboard packaging is used, Peak recommends double boxing of the analyzer using at least 2" of resilient packing material between the inner and outer box walls..

4.0 Peak Performer Operator Interface

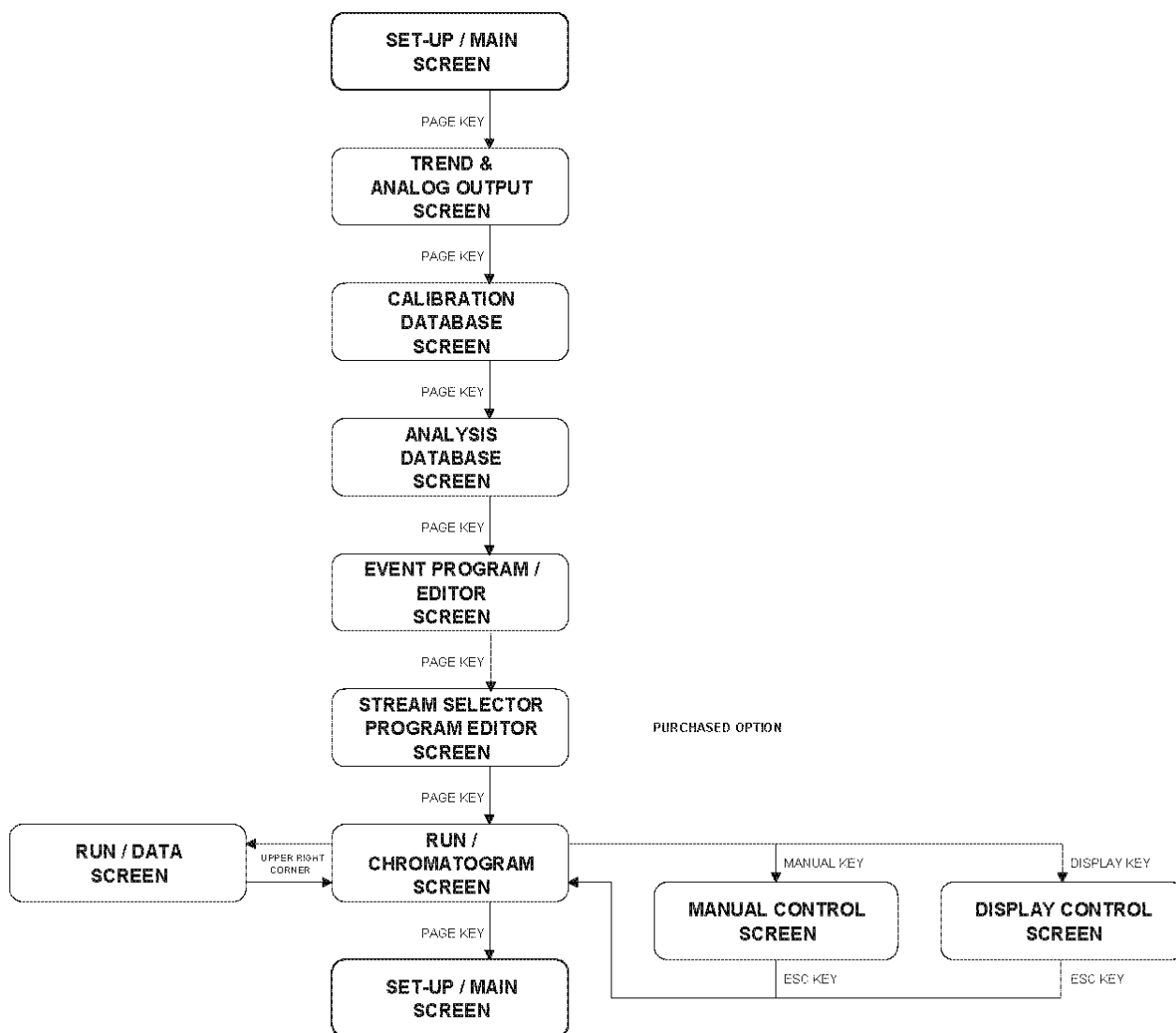


Figure 10. Software Menu Structure

4.1 Set-up Screen

The Set-up screen is multi-functional. Different command buttons appear at the bottom line of the screen dependent upon the section containing the asterick cursor ("*").

1. The upper section (3 lines) contains commands related to analyzer set-up.
2. The middle section (5 lines) contains commands related to temperatures zones and detector parameters.
3. The lower section is reserved for error messages and non-volatile memory updates.

4.1.1 Set-up Screen Upper Section

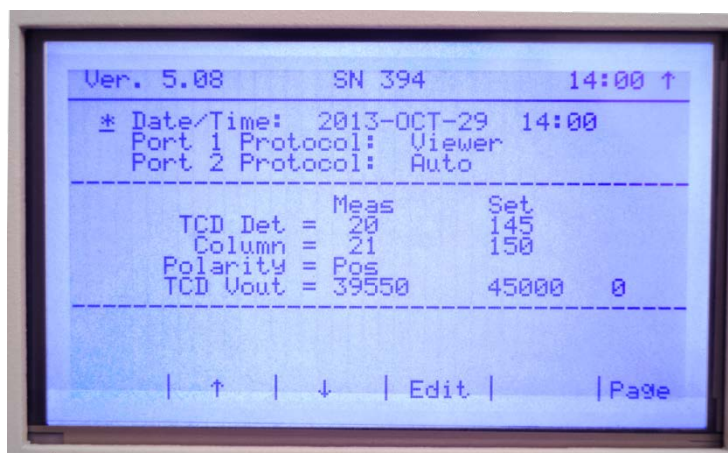


Figure 11. Set-up (or Main) Screen – Upper Section

4.1.1.1 Set-up Screen Upper Section Button Functions

- ↑ : Moves the cursor "*" upwards.
- ↓ : Moves the cursor "*" downwards.
- Edit:** Allows the user to change the Date, Time, or Protocols.
- Page:** Advances the user to the next screen.

4.1.2 Set-up Screen Middle Section (Power-up)

The middle section of the Setup screen is concerns temperature zone and detector function.

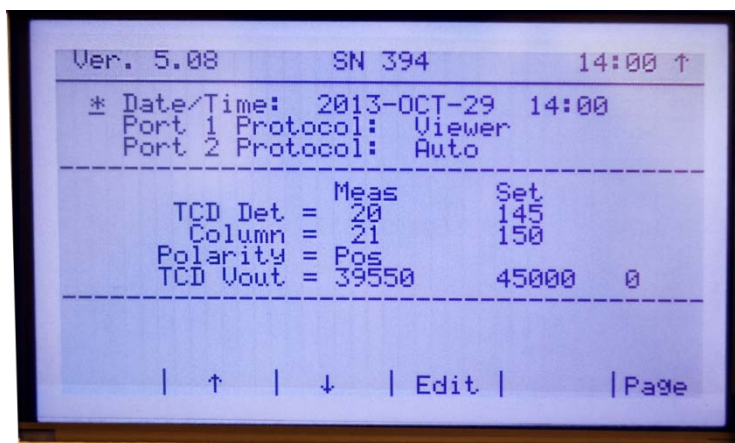


Figure 12. Screen Commands Upon Power-up

NOTE: The heater zones are not energized automatically upon power-up, and the detector electronics are not fully initialized until the ZERO command is invoked.

4.1.2.1 Set-up Screen Middle Section Button Functions (Power-up)

If both heater zones are de-energized, the following command buttons are available:

- | | |
|--------------|--|
| Zero: | Allows the initialize the detector electronics and establish a consistent baseline signal. |
| ↑ : | Moves the cursor "*" upwards. |
| ↓ : | Moves the cursor "*" downwards. |
| Edit: | Allows the user to energize each heater zone, change heater setpoint temperatures, set TCD polarity and set the TCD Vout target value. |
| Heat: | Energizes all the heater zones. |
| Page: | Advances the user to the next screen. |

4.1.3 Set-up Screen Middle Section (Normal)

The middle section of the Setup screen is concerns temperature zone and detector function.

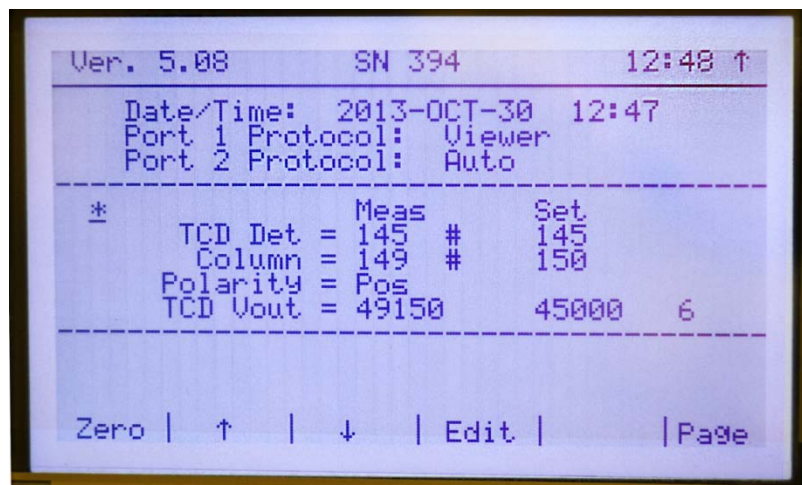


Figure 13. Screen Commands After Warm-up

4.1.3.1 Set-up Screen Middle Section Button Functions (Normal)

Zero: Allows the initialize the detector electronics and establish a consistent baseline signal.

↑: Moves the cursor "*" upwards.

↓: Moves the cursor "*" downwards.

Edit: Allows the user to energize / de-energize each heater zone or change heater set point temperatures, set TCD polarity and set the TCD Vout target value.


Page: Advances the user to the next screen.

4.1.4 Set-up Screen Lower Section

The lower portion of the Setup screen is reserved for error messages and non-volatile memory updates.

4.1.4.1 Error Messages on the Set-up Screen Lower Section

All error conditions are enunciated by a RED color on front panel LED, and a shift in the ERROR FLAG parameter transmitted through the COM1 and COM2 output ports – see Section 6.1 for additional details.

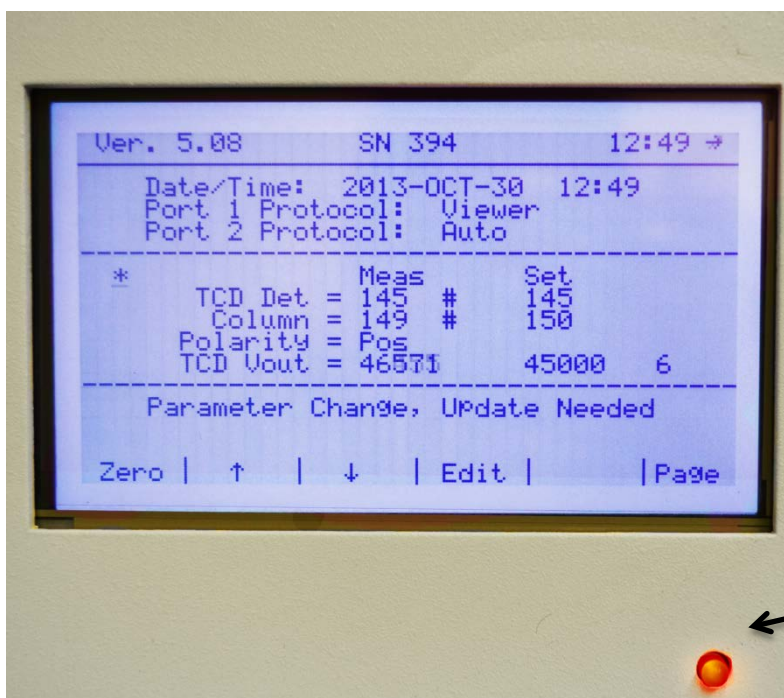

Set-up Screen Displaying an Error Message

Figure 14.

There are eleven individual error conditions that can cause the status LED on the front panel to glow red:

- Detector Communication Error
- Temperature Zone - Out of Range
- Temperature Zone Disabled
- Detector Zeroing Target Error
- Vlamp Low Voltage - Out of Range (RCP analyzers Only)
- FID Flame Temperature - Out of Range (FID Analyzers Only)
- Event Program Load Error
- Stream Selector Program Error
- System Parameters Corrupted
- Parameter Change, Update Needed
- Detector Power Disabled (For PDHID & TCD Analyzers Only)

4.1.4.2 Clearing Error Messages and on the Set-up Screen Lower Section

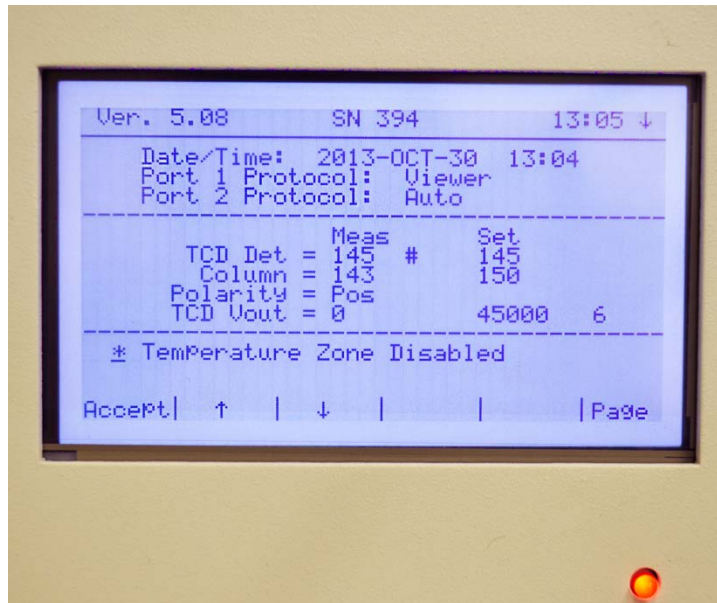


Figure 15. Typical TCD Error Message

All error messages (except PARAMETER UPDATE) can be cleared at any time by pressing the ACCEPT button. However, if the error condition continues to exist, the error message may re-appear.

Accept: Allows the user to acknowledge and clear error messages.

↑ : Moves the cursor "*" upwards.

↓ : Moves the cursor "*" downwards.

Page: Advances the user to the next screen.

4.1.4.2.1 Set-up Screen Lower Section Button Functions

The PARAMETER CHANGE, UPDATE NEEDED error message can only be cleared when the analyzer is not performing an analysis, i.e. is in IDLE mode.

The parameter update action forces a complete re-write of permanent memory and is irreversible.

Peak Labs recommends reviewing all parameter screens before performing a parameter update.

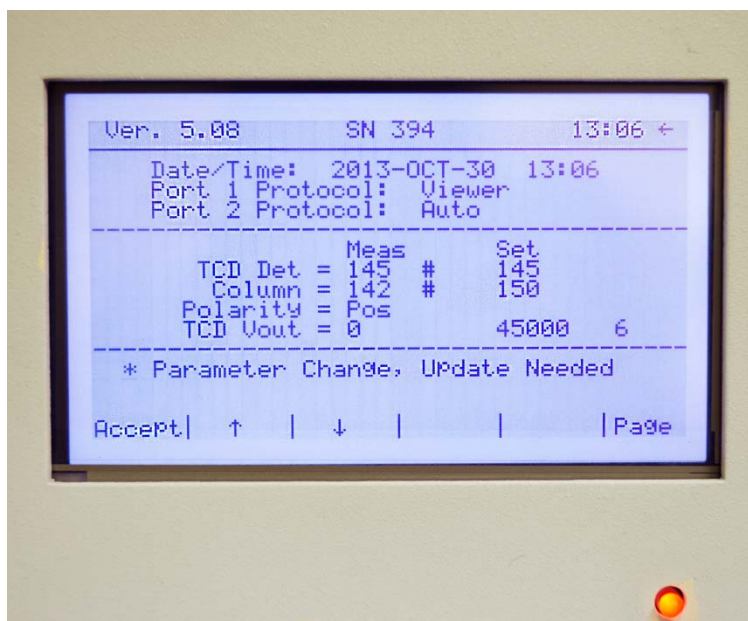


Figure 16. Parameter Change Error Message

Accept: Allows the user to initiate permanent updates to non-volatile memory.

↑: Moves the cursor "*" upwards.

↓: Moves the cursor "*" downwards.

Page: Advances the user to the next screen.

Pressing the ACCEPT button alters the screen appearance.

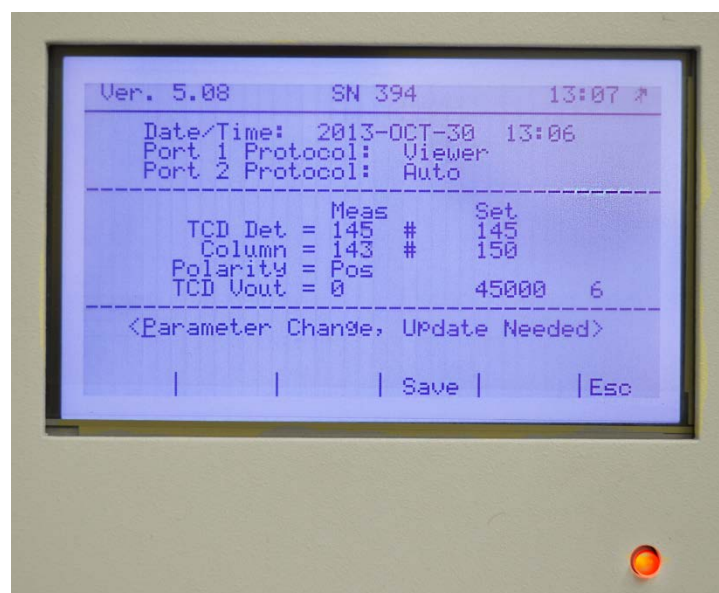


Figure17. Parameter Change Prompt

- Save:** Allows the user to make changes to permanent non-volatile memory.
- Esc:** Aborts updates to non-volatile memory.

4.2 Analog Interface (Trend and Recorder Output) Screen

The compound (or Channel) information for the ANALOG INTERFACE screen is automatically populated using the information from the ANALYSIS DATABASE screen.

Analog signals generated from this information appear on the rear panel's green DIN terminal blocks. See Figure 19.

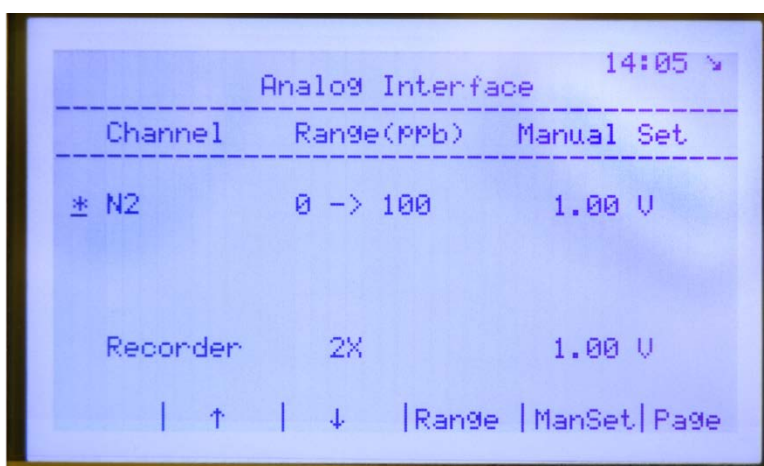


Figure 18. Analog Interface Screen

4.2.1 Analog Interface Screen Button Function

↑ : Moves the cursor "*" upwards.

↓ : Moves the cursor "*" downwards.

Range: Allows the user to select the full scale concentration range associated with the maximum 1.00 VDC signal for each compound. The ranges are manually set from 0 ->999999 ppb

e.g. 0 – 1000 ppb readings in the analyzer = 0 – 1.00 VDC output scale, where a 500 ppb value would output 0.5 VDC.

At the end of each run, the analog signal is updated and holds until the next run is completed. As an accessory, Peak can provide standard commercial VDC -> mA converters as needed.

In lieu of a specific concentrations range, the RANGE button for the recorder function toggles through choices of raw signal attenuation: X1, X2, X4, X8, X16, X64. This function is similar to the RANGE switch on a chart recorder. This function only affects the recorder output not the individual compound trends.

Manset: Temporarily over-rides the signal output with a known output voltage for test purposes.

Four MANSET options are available: 0.00 VDC, 0.25 VDC, 0.5 VDC, and 1 VDC. There is not any zero offset or full voltage scale adjustment within the PP1 – if this function is a requirement, Peak Labs recommends the use of standard commercial voltage converters attached externally to the analyzer.

Page: Advances the user to the next screen.

4.2.2 Analog Interface Connections

All analog interface connections are made via the rear panel DIN connector. To make a proper connection, strip all wires back approx. 3/8", unscrew the terminal locking screw fully, insert the wire and tighten. A common error is insufficient length of bare wire creating the condition where the terminal contacts only the wire insulation.



ANALOG 1	⊙	⊙	COMMON GND
ANALOG 2	⊙	⊙	COMMON GND
ANALOG 3	⊙	⊙	COMMON GND
ANALOG 4	⊙	⊙	COMMON GND
ANALOG 5	⊙	⊙	COMMON GND
ANALOG 6	⊙	⊙	COMMON GND
RECORDER	⊙	⊙	COMMON GND
REMOTE START	⊙	⊙	COMMON GND
FID FLAME +	⊙	⊙	FID FLAME -
GEN. ALARM +	⊙	⊙	GEN. ALARM -

Figure 19. Analog Output Connections



NOTE: The analog outputs are locked if the Port 1 Protocol on the MAIN SCREEN is set to DETECTOR. See Section 6.2.

4.2.3 Analog Interface Wiring Breakdown

There are 10 total terminal connectors, from the top, the connectors ID's are:

Terminals 1- 6:

- Impurity Trends:
All channels are mapped to the calibration database
All channels are 0 – 1 VDC
All channels signals are updated at the end of run and held
All channels are manually scalable on the "ANALOG INTERFACE SCREEN"
e.g. 0 → 1000 = 0 to 1000 ppb = 0 – 1 VDC scale. Maximum scale is 999,999
Optional converters to provide other current / voltage outputs are available. Please contact your local Peak representative for additional information.

Terminal 7:

- Analog recorder / raw signal
0 – 1 VDC range

Terminal 8:

- Remote Start Input:
The remote start input is activated by shorting the terminal 8 to COMMON GND for at least 100ms. The remote start input must then return to an open state (unshorten).

Terminal 9:

- FID Flame Switch Relay:
The FID Flame and General Alarm switch relays are SPST- N.O. and rated for 0.5 amps @ 200VDC (10VA)

Terminal 10

- General Alarm Relay:
The FID Flame and General Alarm switch relays are SPST-N.O. and rated 0.5 amps @ 200VDC (10VA)



NOTE: The Analog recorder signal is only active when the COM1 setting is set to "Viewer".


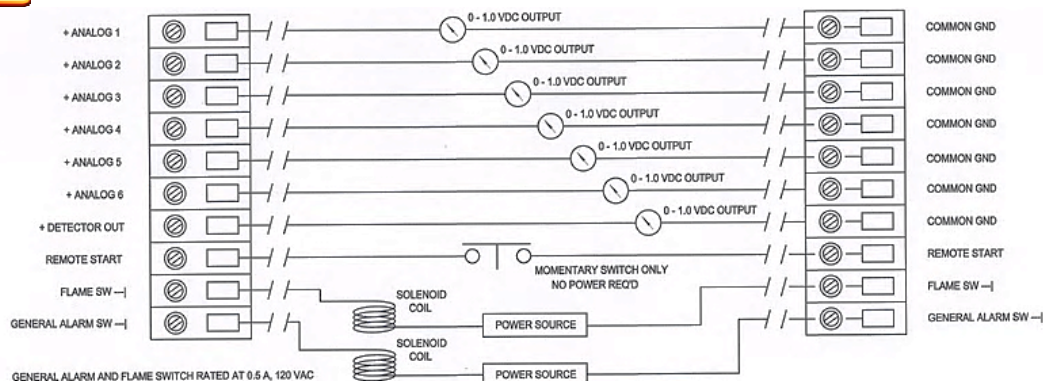


Figure 20. Analog Output Wiring Schematic

Channels 1-7 output voltage can be tested manually by pressing the MANSET key to send a temporary signal for diagnostics.

4.3 Calibration Database Screen

The Calibration Database screen stores critical information regarding calibration (span) gas concentrations and instrument response, and provides tools to easily update the response factors.

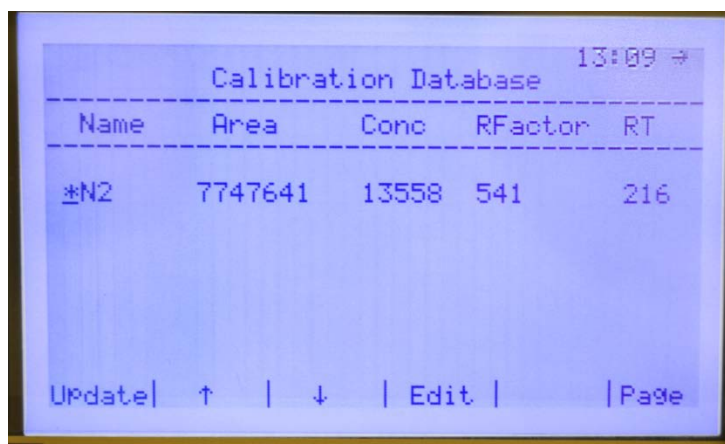

Calibration Database 13:09 →				
Name	Area	Conc	RFactor	RT
*N2	7747641	13558	541	216

Update | ↑ | ↓ | Edit | Page

Figure 21. Calibration Database Screen

4.3.1 Calibration Screen Button Function

Update: Allows the analyzer to calculate new RF Factor.

↑ : Moves the cursor "*" upwards.

↓ : Moves the cursor "*" downwards.

Edit: Allows the user to change concentration data or response factor.

Page: Advances the user to the next screen.

4.3.2 Calibration

The Peak Performer 1 chromatograph determines concentration peak area to quantify the amount of contaminants in the sample gas. The direct instrument response (in units of uV-sec) is reported simultaneously as an overlay on the RUN / CHROMATOGRAM SCREEN and in the AREA column on the CALIBRATION DATABASE SCREEN.

As the concentration peak area is determined, the analyzer applies the pre-programmed response factor as an inverse proportional factor to generate the reported concentration value.

$$\text{PEAK CONCENTRATION} = \text{PEAK AREA RESPONSE} / \text{RESPONSE FACTOR}$$

4.3.2.1 Establishing New Instrument Response Factors

Connect a known calibration standard gas (span gas) to the analyzer at the SAMPLE IN port on the front panel. (If span gas connection is made at the rear panel SAMPLE IN port, ensure that the front jumper is connected between SAMPLE OUT and SAMPLE IN on the front panel)

Verify that sample is flowing freely through loop - exit flow from the SAMPLE OUT should be 30 to 120 cc/min.

Press the ESC or PAGE keys as needed to navigate to the RUN screen

Press the SINGLE / CYCLE / RERUN key to toggle to SINGLE mode

Press the RUN key to initiate analysis – the P1 event program will be loaded automatically and executed. Allow the analyzer to fully execute the event program. Current concentration data will be displayed as an overlay on the RUN screen as peaks are identified and quantified.

Press the PAGE key as needed to toggle to the CALIBRATION DATA BASE screen

Place the “ * ” at the first compound of interest and press the EDIT key – use the up, down, left and right arrow keys to enter the span gas concentration in the CONC field. Be certain to use the right arrow key to fully exit the CONC field or all updates will be lost.

With the “ * ” before the compound name, press the UPDATE key. The analyzer will recalculate the appropriate response factor for the compound based upon the span gas concentration and analyzer area response using the following formula:

$$\text{RESPONSE FACTOR} = \text{AREA RESPONSE} / \text{SPAN GAS CONCENTRATION}$$


NOTE: The LED on the analyzer's front panel will change from green to red to remind the user that the current calibration information is stored in volatile (temporary) memory only, and will need to be saved to permanent memory at a later time.

Repeat the sequence for each of the compounds shown on the CALIBRATION DATA BASE screen.

Repeat analysis of the span gas sample and verify that the analyzer reports data values within normal calibration tolerances (i.e. + / - 5%).

4.4 Analysis Database Screen

The Analysis Database screen stores critical information regarding peak identification and quantification parameters.

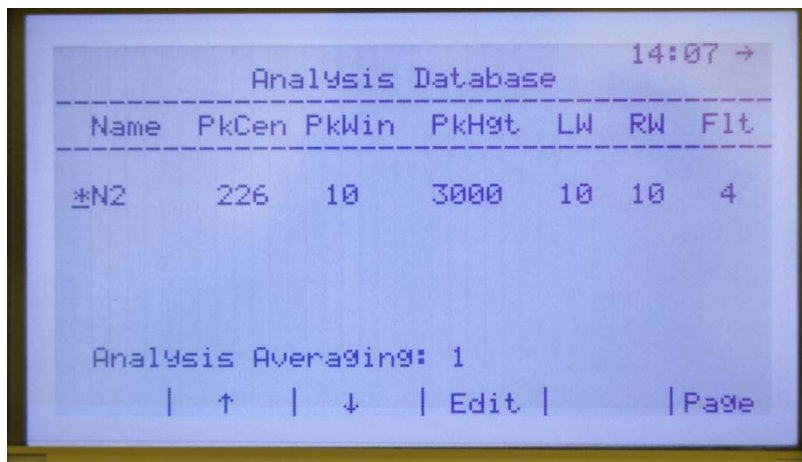

The screenshot shows a terminal window titled "Analysis Database" with a timestamp "14:07" in the top right. The table contains one data row for peak "*N2". Below the table, it shows "Analysis Averaging: 1" and a control bar with navigation icons and labels.

Name	PkCen	PkWin	PkHgt	LW	RW	Flt
*N2	226	10	3000	10	10	4

Analysis Averaging: 1

| ↑ | ↓ | Edit | Page

Figure 22. Analysis Database Screen

4.4.1 Analysis Database Parameters

Name = The peak identification tag, expressed as 5 alphanumeric characters (A-Z, 0-9)

PkCen = The normally expected peak retention time (in seconds), used for assigning the appropriate "Name"

LW = The typical time span measured from the start of the peak's baseline rise to the peak apex "PkCen" (in seconds). Minimum value is 3, maximum is 35

RW = The typical time span measured from the peak's apex "PkCen" to the end of the peak's baseline decline (in seconds) Minimum value is 3, maximum is 35

PkWin = Total tolerance window (in seconds) for assignment of a "Name" to a quantified chromatographic peak. The tolerance window is centered upon the "PkCen" value

PkHgt = A variable for establishing the cross-over point in peak detection modes. Chromatographic peaks higher than this parameter will be quantified using the "Variable" mode heights, peaks with height less than or equal to "PkHgt" will be quantified using the "ForceB" mode.

Flt = The convolution filter value expressing the overall peak shape. Flt = 2 is recommended for sharp, narrow peaks such as H₂, and Flt = 8 is recommended for broad peaks such as CO.

As a general rule, LW and RW values do not change significantly with peak concentration, so calibration peak size is not critical in determining proper LW and RW values.

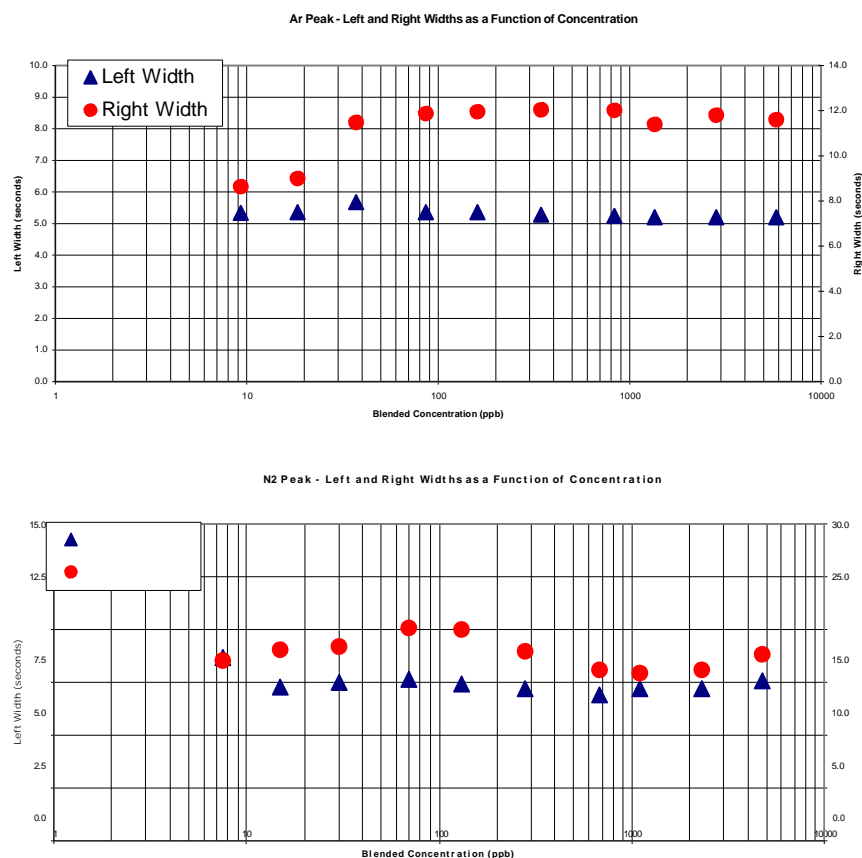


Figure 23. Left and Right Peak Width as a Function of Concentration

- Total peak width (LW + RW) is primarily a function of sample loop unloading.
- Larger sample loops require a long time period to fully inject.
- The relative values of LW and RW are dictated by the skew factor of the peak.
- Peaks need about 25 –30 data points minimum to be well measured - at 5 data points per second, this means total peak width should be greater than 6 seconds.

Therefore, it is recommended to use a substantial calibration peak concentration to facilitate easy LW and RW measurements.

4.4.2 Analysis Screen Button Function

↑ : Moves the cursor “*” upwards.

↓ : Moves the cursor “*” downwards.

Edit: Allows the user to enter new peak compound name or analysis parameters.

Page: Advances the user to the next screen

Analysis Averaging - When enabled, one, two or four chromatograms are averaged (data point by data point) to create a composite chromatogram, which is then quantified using normal techniques. Analysis averaging does not average concentration data values.

From analysis start, it will take approximately 30 minutes to obtain the best averaged results.

4.5 Event Program Editor Screen

The Event Program Editor screen stores critical information regarding pre-programmed instrument actions needed to perform an analysis.

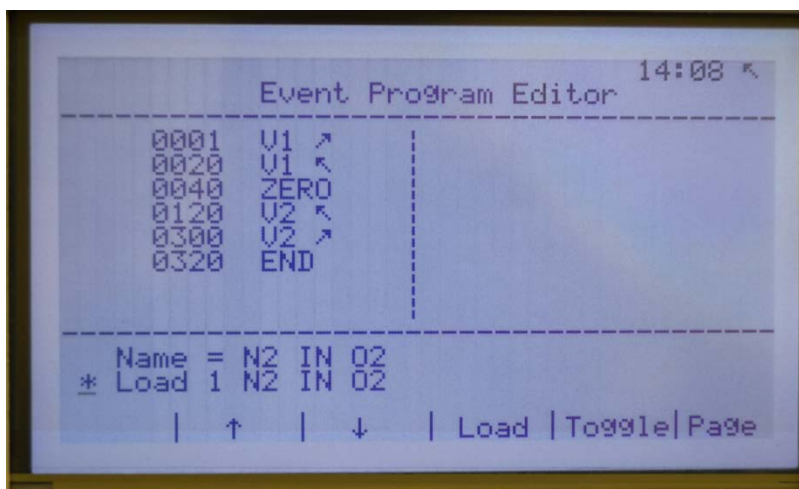


Figure 24: Event Editor Screen

4.5.1 Event Program Commands

The instrument actions can be pre-programmed to create an analysis method or Event Program. The analyzer stores up to 4 distinct Event Programs.

V1/ (CW)	Moves V1 to the Inject position
V1 \ (CCW)	Moves V1 to Load position
V2/ (CW)	Moves V2 to the In-line position
V2 \ (CCW)	Moves V2 to the Vent position
ZERO	Invokes the detector ZERO command
END	Ends the Event Program and prepares the analyzer for the next run
P ON	Not currently in use
P OFF	Not currently in use

4.5.2 Event Editor Screen Button Function

- ↑ : Moves the cursor "*" upwards.
- ↓ : Moves the cursor "*" downwards.
- Load 1:** Allows the user to load the event program # 1.
- Toggle:** Allows the user to scroll through event program choices.
- Page:** Advances the user to the next screen

4.6 Stream Selector Program Editor Screen (Purchased Option)

Peak offers a dual sample stream option for most analyzers. The selection and timing of sample streams can be pre-programmed on this screen.

NOTE: This screen will only appear when the STREAM method of analysis is selected on the RUN screen

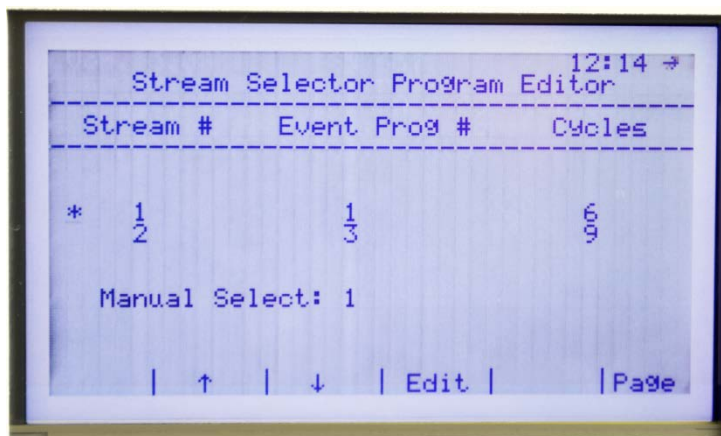


Figure 25. Stream Selector Sequence Screen

4.6.1 Stream Selector Commands

The selection of sample streams and appropriate Event Program can be run in a pre-programmed sequence.

- | | |
|--------------|--|
| Stream # | Establishes which sample stream to be utilized for analysis |
| Event Prog # | Establishes which Event Program to be utilized for the sample stream |
| Cycles | Establishes the number of repeat analysis before switching to other stream |

4.6.2 Stream Selector Program Editor Screen Button Function

- ↑ : Moves the cursor "*" upwards.
- ↓ : Moves the cursor "*" downwards.
- Edit:** Allows the user to chose event program # and # of repetitions.
- Page:** Advances the user to the next screen

4.7 Run / Chromatogram Screen

The Run / Chromatogram screen provides the real time display and data processing of sample analysis. The screen displays the running chromatogram, peak markers, and concentration data. The information displayed on this screen is useful for evaluating instrument health, event program validity, and analytical stability,

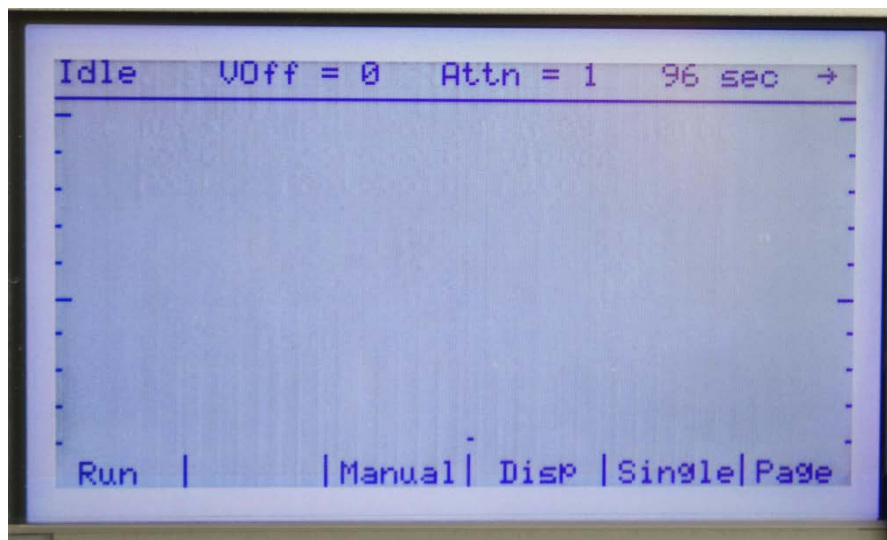


Figure 26. Run / Chromatogram screen in the idle state, and the mode displayed upon initial power up

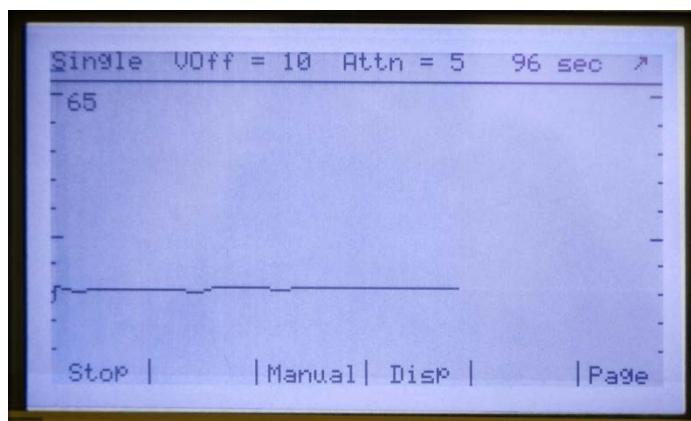


Figure 27. Run screen in Single mode, with chromatogram being drawn

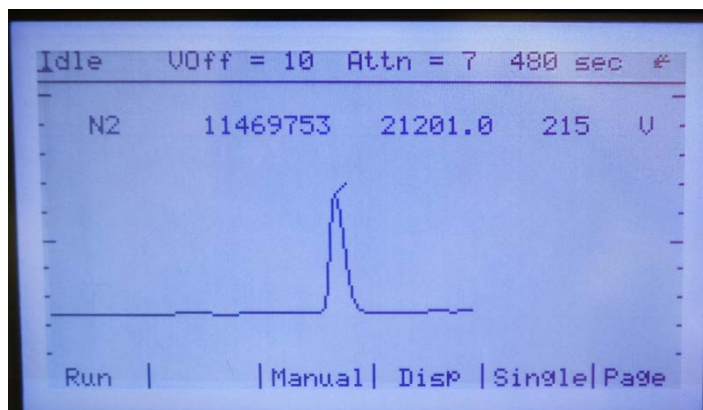


Figure 28. Run screen in Idle mode at end of run

4.7.1 Run Screen Button Functions

- Run:** Loads a standard event program and runs an analysis in the mode selected.
- Idle:** Stops the analysis immediately
- Manual:** Opens the Manual Run Screen
- Disp:** Opens the Display Screen
- Single:** (Toggle). Starting a run in the Single mode will command the PP1 to run one time and return to the idle mode.

- Cycle:** (Toggle). Starting a run in the Cycle mode will configure the PP1 to run continuously until the Idle key is pressed.
- ReRun:** (Toggle). Starting a run in the ReRun mode will recompute the area and concentration based on the parameters in the analysis data base.
- Stream:** (Toggle). Starting a run in the Stream mode will configure the PP1 to initiate the Stream Selector Sequence until the Idle key is pressed.
- Page:** Advances the user to the next screen.

4.7.2 Continuous Monitoring Analysis

Continuous monitoring is intended for on-line analysis or other applications which ReRun analysis is not required.

Use the "PAGE" key to scroll to reach the Run screen. Toggle the Single / Cycle / ReRun key until Cycle appears. (see Figure 28) Press Run to initiate a repetitive analytical cycle

By pressing the "RUN" key, event program 1 (P1) is automatically loaded and the system will begin to collect data in repetitive mode (see Figure 28).

To halt the cycle, pressing the "IDLE" will cause the analysis run to cease immediately. Toggle between the Single, Cycle, Rerun (or optional Stream) to select the next mode of analysis.

4.7.3 Manual Screen (Sub-Screen of the Run Screen)

The Manual sub-screen allows the operator to invoke all the event program commands on demand.

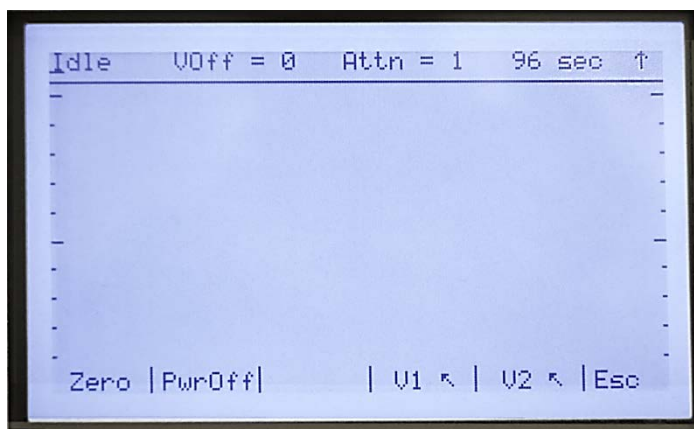


Figure 29. Run / Manual screen

4.7.3.1 Manual Screen Button Functions

Zero: Pressing this key will electronically zero the detector.

PwrOff/On: Energizes the TCD filament current.

V1: Toggles Valve 1 between CW and CCW positions.

V2: Toggles Valve 2 between CW and CCW positions.

Esc: Exits from the Manual sub-screen to the Run screen.

4.7.4 Display Screen (Sub-Screen of the Run Screen)

The Display sub-screen allows the operator to adjust the horizontal scale, vertical scale and baseline offset of the chromatogram. These functions mimic those of a standard chart recorder.

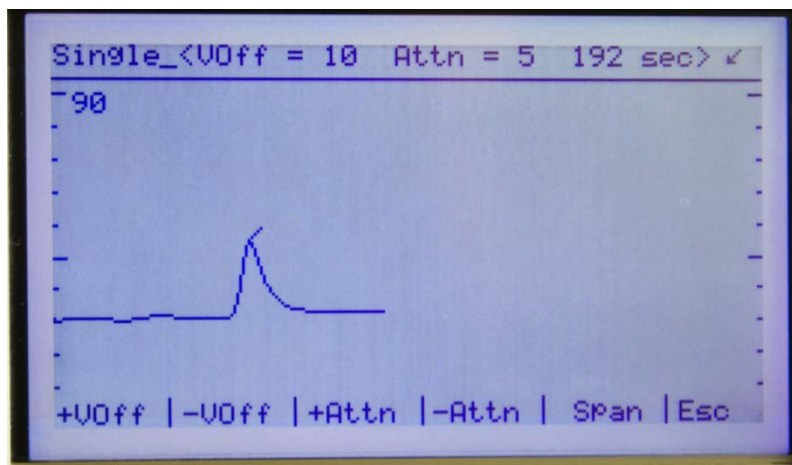


Figure 30. Run / Display screen

4.7.4.1 Display Screen Button Functions

+VOff: Scrolls the screen up (maximum Voff = 11) without changing the scale.

-VOff: Scrolls the screen down (minimum Voff = 0) without changing the scale.

+Attn: Toggles the attenuation / vertical scaling (max. Attn = 11) up.

-Attn: Toggles the attenuation / vertical scaling (min. Attn = 1) down.

Span: Toggles the window view (horizontal scaling) in seconds.

Esc: Exits from the Display sub-screen to the Run screen.

5.0 Peak View Software

The Peak View software has two main capabilities. Primary functions concern acquisition and archiving of compound concentration and acquisition and collection of the chromatogram data.

The software also provides the ability to display and print the data from a remote computer platform.

Peak View software is written specifically for the Peak Laboratories Peak Performer 1 (PP1) gas analyzer covering communications port protocols and data formats and is not intended for use with other manufacturer's analyzers.

5.1 System Requirements

- 1.2 GHz, 20 GB, 256 Mb RAM, 1 USB port
- Windows XP or Windows Vista, Windows 7, Windows 8
- USB to DB9 (RS232) Adapter – Maker: Gigaware Model # 26-949 or equivalent
- DB9 Extension cable, DB9M to DB9F, wired straight-through, 10 feet
Optional : USB 4 port self-powered Hub – Inland # 08302 or equivalent

5.2 Installing the Peak View Software

- 1) Place Peak View software CD-ROM disk into CD-ROM drive.
- 2) If Autorun is enabled, follow the installation screen prompts.
- 3) If the CD-ROM does not automatically launch, select "SETUP" from CD-ROM disk.
- 4) The installation software will auto-install the Viewer in a folder under the PC's main drive as the default. The user has options from the installation process to create a new folder and install in a new location.

5.3 Starting the Peak View Software

- 1) Double-click the Peak View icon on the Desktop or
- 2) Select "START" button on the bottom Taskbar.
- 3) Select "PROGRAMS", select "PEAK VIEW ", and double-click to launch move to "Peak View"

5.4 COM1 / COM2 Port and Cable Wiring

The cable connected from the PC to the PP1 should be a straight-through wired 9 pin M/F cable, attached to the PC's COM port and the PP1's COM port. (See Section 6 for details on communications protocols.)

5.5 Viewer Software Window

The Viewer window shown in Figure 33 shows the main Viewer screen **with** a RED DOT in the upper right corner. This RED DOT signifies that there is a connection problem with the PC and the PP1.

Verifying the following can solve this connection problem.

- Communication Settings
- Connection of 9-pin cable between PC and PP1.
- RS-232 cable plugged into the PP1's COM 1 port.
- PP1 is powered off

After the connection problem has found and corrected, the RED DOT should be **CLICKED** to verify if communication between the PC and PP1 can be properly established.

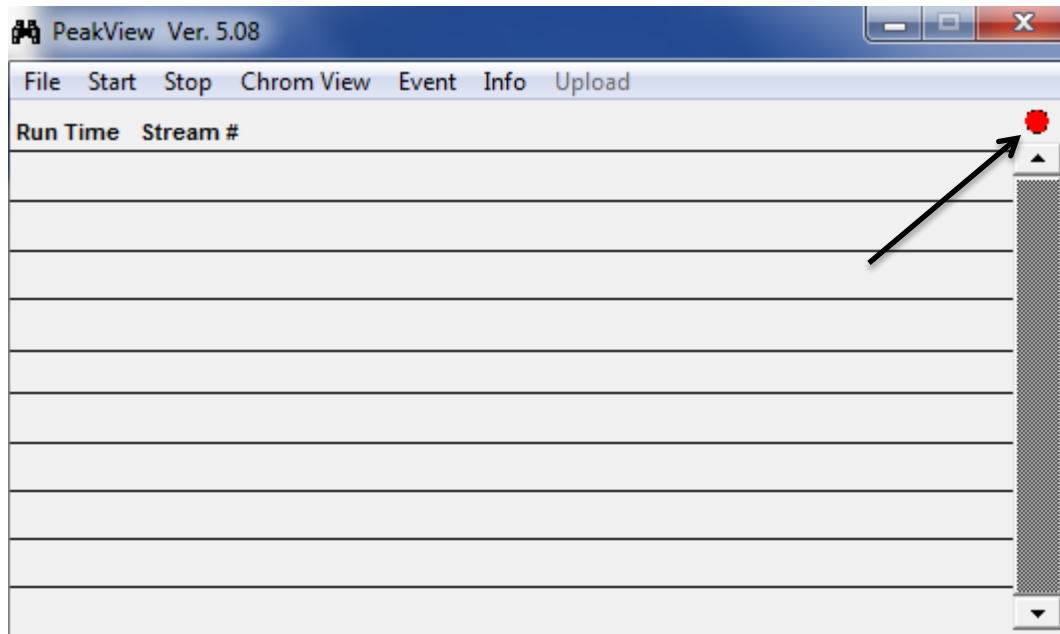


Figure 31. Viewer main window showing Communications Error

The Viewer window shown in Figure 32 shows the viewer screen **without** the RED DOT. This signifies the communication connection between the PC and the PP1 is function correctly.

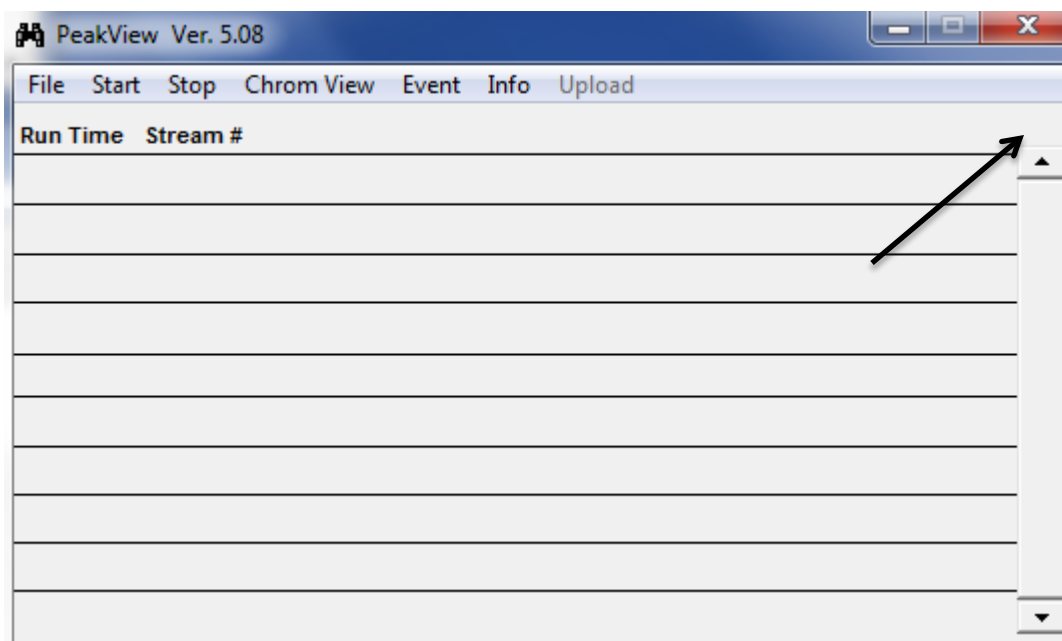


Figure 32. Viewer main window showing no Communications Error

5.5.1 File Menu on the Viewer Taskbar

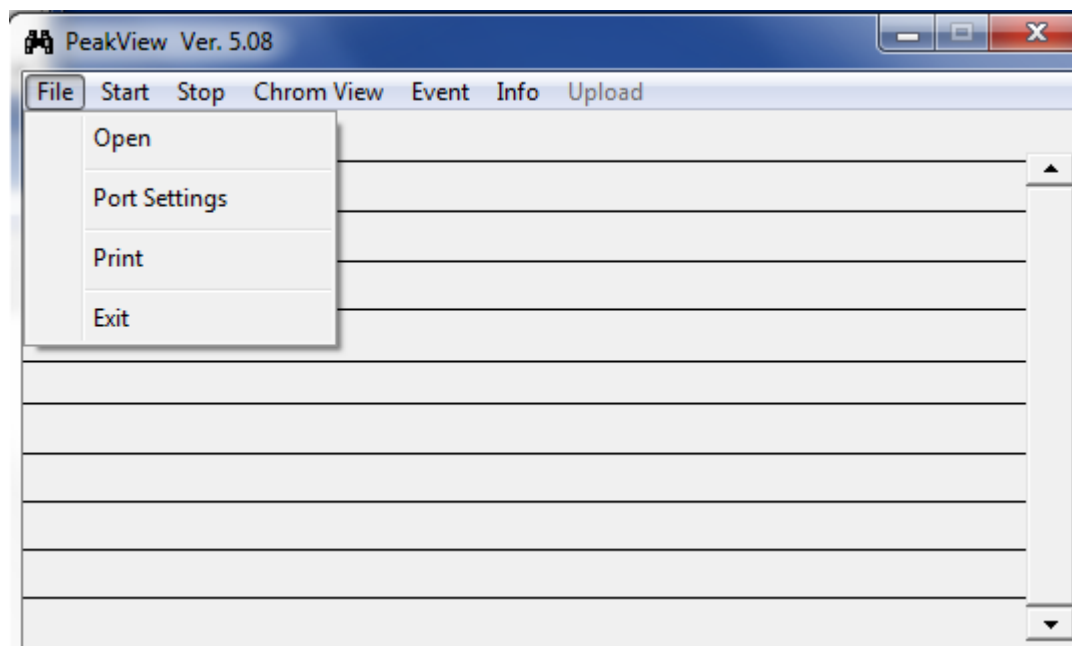


Figure 33. File Menu on the taskbar

The File menu contains the following commands:

Function	Description
Open	Opens a *.CSV file that was created by the Viewer software. Use this command to open a *.CSV file to view the concentration data.
Port Settings	Opens serial communication settings window for setting RS-232 configurations between the Host PC and the PP1. As shown in figure 34.
Print	Prints all concentration data of selected *.CSV file.
Exit	Exit Viewer application.

5.5.1.1 Communication Setting Window

The Communication Settings window is shown in figure 34. This screen displays the serial communication settings used by the PC to communicate to the PP1.

The baud rate, data bits, and stop bits settings should match all of the settings in figure 16, baud rate 38400, data bit 8, stop bit 1.

The COM port number must match the COM port designated by the PC's terminal. The COM Port number is determined by which COM port the serial cable is plugged into on the PC, not the PP1.

Clicking the "OK" button will configure the COM port to the correct options and save the settings in a file on the computer.


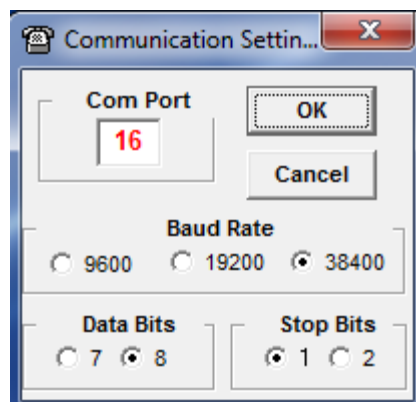


Figure 34. Viewer Software Communications Setting screen

Once the *.CSV file is opened, the Viewer window automatically updates and shows concentration data. Figure 35 is an example with consecutive runs. The third run of concentration data at 08:43 is displayed in RED to denote that there was a general error in the PP1 during this run. The error was cleared before the end of the next run at 08:51.

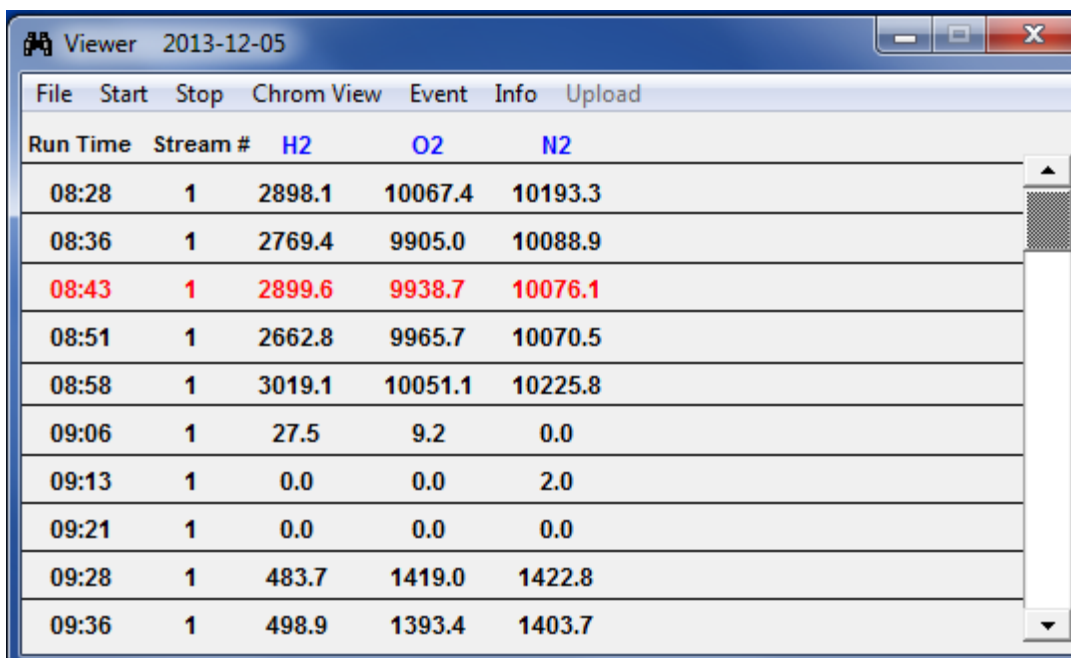

Run Time	Stream #	H2	O2	N2
08:28	1	2898.1	10067.4	10193.3
08:36	1	2769.4	9905.0	10088.9
08:43	1	2899.6	9938.7	10076.1
08:51	1	2662.8	9965.7	10070.5
08:58	1	3019.1	10051.1	10225.8
09:06	1	27.5	9.2	0.0
09:13	1	0.0	0.0	2.0
09:21	1	0.0	0.0	0.0
09:28	1	483.7	1419.0	1422.8
09:36	1	498.9	1393.4	1403.7

Figure 35. Viewer screen showing concentration data

There are eleven individual error conditions that can cause the Viewer Line to appear red:

- Detector Communication Error
- Temperature Zone - Out of Range
- Temperature Zone Disabled
- Detector Zeroing Target Error
- Vlamp Low Voltage - Out of Range (RCP Analyzers Only)
- FID Flame Temperature - Out of Range (FID Analyzers Only)
- Event Program Load Error
- Stream Selector Program Error
- System Parameters Corrupted
- Parameter Change, Update Needed
- Detector Power

5.5.1.2 Viewer Line Details Window

This window is displayed by double clicking on any time-stamp on the Viewer window. The screen provides detailed analytical information about the run selected.

Viewer Details 2013-12-05 08:51

Print

Name	H2	O2	N2			
Peak Time (sec)	74.4	106	140.2			
Left Width (sec)	5.8	8.8	8.8			
Right Width (sec)	7.2	26.8	24			
Measured Area	157103	27774508	29476416			
Response Factor	59	2787	2927			
Width Analysis	Variable	Variable	Variable	N/A	N/A	N/A
Flame Temperature	0.0	Run Type	Cycle	Stream # 1		
System Alarm	OK	Analysis Index	3			

Figure 36. Viewer Detail screen showing raw analysis data.
(The information on the screen is linked to run 08:43 displayed in figure 35.)

5.5.2 Start Menu on the Viewer Taskbar

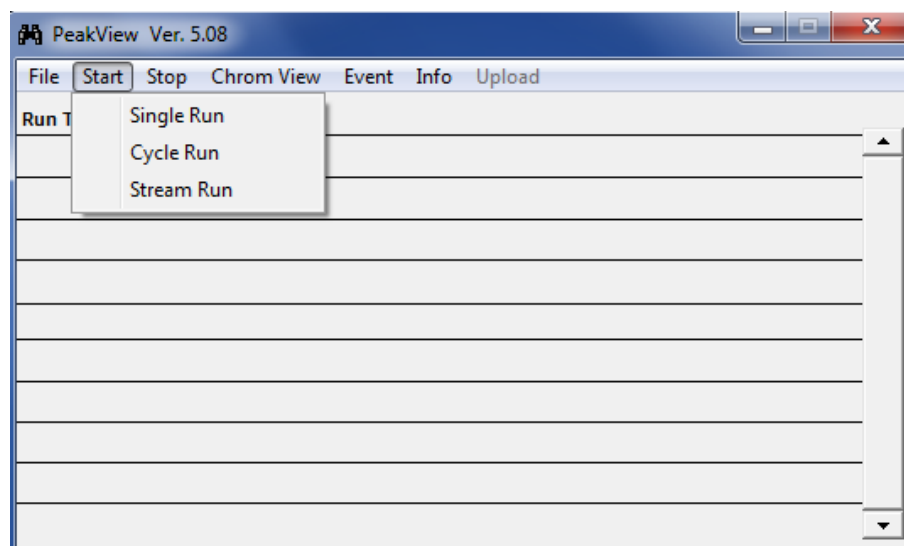


Figure 37 Start menu on the taskbar

The Start menu contains the following commands:

Function	Description
Single Run	Single Run selection commands the PP1 to start a Single Run.
Cycle Run	Cycle Run option commands the PP1 to start a continuous Cycle Run.
Stream Run	Stream Run selection commands the PP1 to start a Stream Selection Run. (This operation only function if the Stream Selection option is installed)

5.5.3 Stop Menu on the Viewer Taskbar

The Stop Command selection instantly sends a command the PP1 to terminate the current run immediately.

5.5.4 Chrom View Menu and Chromatograph Viewer Window

The Chrom View menu selection opens the Chromatograph Viewer window. Two methods can open the Chrom Viewer window.

- 1) The first method is to select the Chrom View menu selection from the Viewer taskbar and then select a ".chm" file to display.
- 2) The second method of opening the Chrom Viewer is double-clicking one of the concentration values displayed in the Viewer window.

Using either method, after the chromatogram data is loaded, the display can be adjusted by the Offset, Attenuation, and Span scroll bars.

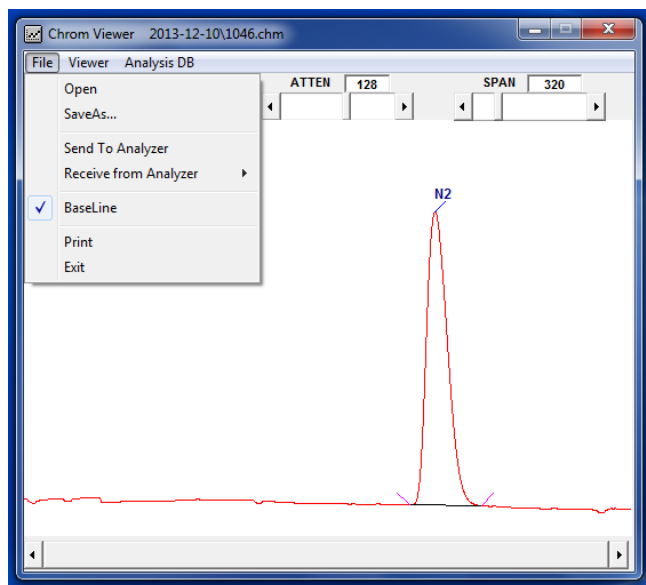


Figure 38. Chromatogram Viewer Window

NOTE: Viewer window scaling does not match the analyzer's display scaling . A higher degree of resolution is available in the Viewer window. For instance, Viewer ATTN 8 approximately is equivalent to analyzer display ATTN 1

5.5.4.1 Chromatogram Files

Chromatograms are archived on the host computer as text files, suffixed ".chm". The text files can be imported into most common laboratory analysis, spreadsheet and word processor programs as needed.

Data points in the chromatogram file are time sequenced at 0.2 seconds per data point.

The filename of the chromatogram file is generated by time of day on the host computer, and the file is stored in a folder that corresponds to the date of collection.

A corresponding Viewer data file is created and / or updated at the time the chromatogram file is created.

As an example, the fourth chromatogram shown in Figure 36 would be archived on the host computer as follows:

Main directory	C:\ Viewer
Sub-directory	2013-12-05 [(YYYY-MM-DD) of archive creation]
Viewer filename	2013-12-05.csv
Filename	0851.chm

5.5.4.2 File Menu on the Chrom Viewer Taskbar

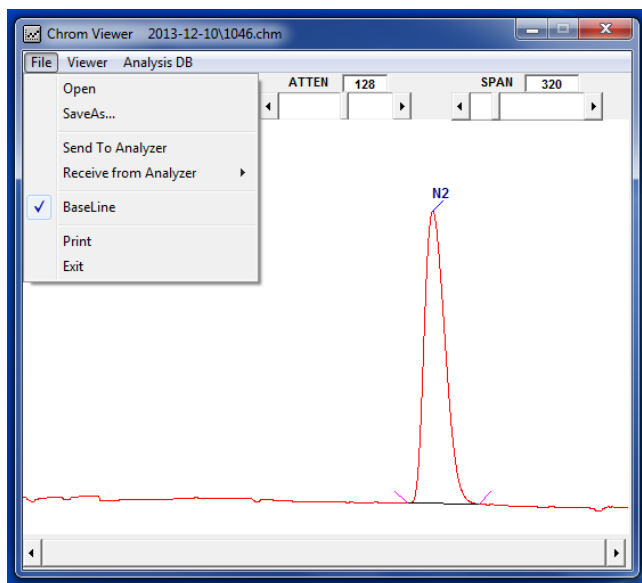


Figure 39. Baseline Selected

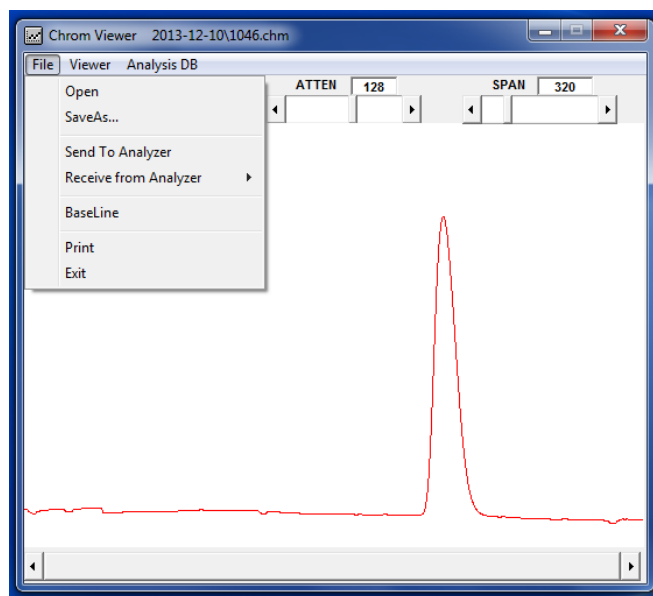


Figure 40. Baseline Deselected

The File menu offers the following commands:

Function	Description
Open	Opens a chromatogram file. Use this command to open and load the chromatogram data for display.
Save As	Save the chromatogram data that is currently displayed.
BaseLine	Enables display of compound name and integration baseline information.
Print	Prints window.
Exit	Exit window.

5.5.4.3 Viewer Menu on the Chrom Viewer Taskbar


The Viewer menu selection returns the user to the original Viewer window. (See Figure 36)

6.0 Communication Port Protocol and Formats

There are two serial communications ports on the rear of the PP1. COM1 and COM2 are multifunctional serial ports.

6.1 COM1 / COM2 Port and Cable Wiring

The cable connected from the PC to the PP1 should be a straight-through wired 9 pin M/F cable, attached to the PC's COM port and the PP1's COM port. Figure 42 is the minimum cable configuration.

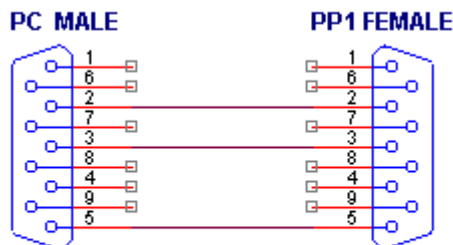

9 Pin Cable from PC to PP1

Figure 41. COM Cable Wiring

6.2 Port 1 Protocol

There are two COM1 output formats: Peak View and Detector. Proper selection of protocol is determined by the application in use on the host computer.

6.2.1 COM 1 Port Settings

The RS-232 settings for COM 1 port are 38400 baud rate, 8 bit data, 1 stop bit, and no parity. These settings are fixed and are not user-selectable.

6.2.2 Peak View Output (Columns Format in .CSV files)

The Peak View software stores information in a common text format, Comma Separated Variable (CSV). CSV files can be opened in common spreadsheet applications for further analysis.

The column format in the *.CSV files are:

A = Data/ Time

B = Run Mode

C = Analysis Buffer Number and Stream Number

D = General Error Alarm

E = RCP Lamp Voltage or FID Flame Temperature

(NOTE: No values to be measured for TCD or PDD and will show up as 0)

F = Compound Name

G = Retention Time

H = Left Point Offset

I = Right Point Offset

J = Area

K = Response Factor

L = Processing Flag- V,F or B

Columns F to L are repeated for M to S, T to Z, and AA to AG for the next five compounds listed in the Calibration Data Base.

6.2.3 Detector Output

The detector output of the COM1 interface currently supports **one** command from a host computer system. This command remotely starts a RUN.

The PP1 system processes the external start command only if the analyzer is in IDLE mode. (Not running a sample). Upon receipt of the external start command, the PP1 performs the following functions:

- Sets system run mode to SINGLE RUN mode.
- Loads event program P1.
- Starts a Single Run.
- Sends Detector Signal Data stream (i.e. Vout reading) to the COM1 port.
- At the end of the run, the PP1 sends an End-Of-Data string to identify the end of the data stream.
- Sets analyzer to IDLE mode and waits for the next start command.

6.2.3.1 External Start Command String

The External Start Command string format transmitted by host computer must consist of three ASCII characters. Below are four forms of the same command:

Form # 1	Start-Of-Text "S" End-Of-Text
Form # 2	STX "S" EXT
Form # 3	^B S ^C
Form # 4	02 53 03

6.2.3.2 Detector Signal Data String

The PP1 Detector Signal Data string format is an integer number followed by a carriage return (CR) and line feed (LF).

An example:

```
104345
104763
104832
...
```

The End-Of-Data string format consists of three ASCII characters. The following are four forms of the same command:

Form # 1	Start-Of-Text "E" End-Of-Text
Form # 2	STX "E" EXT
Form # 3	^B E ^C
Form # 4	02 45 03

An example:

```
...
104345
104763
104832    // Last sample sent to Host
^BE^C    // End-Of-Data command to Host from PP1
```

6.3 Port 2 Protocol

There are two COM2 protocol formats, PLC (programmable logic controller) and MODBUS protocols. There are three modes of operation: AUTO, POLL, and MODBUS. PLC protocol can operate in 2 modes, AUTO or POLL. The MODBUS Protocol mode is predefined and not user changeable. Proper selection of protocol is determined by the application in use on the host computer.

6.3.1 COM 2 Settings

The RS-232 settings for COM 2 port are 9600 baud rate, 8 bit data, 1 stop bit, and no parity for AUTO and POLL modes. The RS-232 settings for COM 2 port in MODBUS mode are 9600 baud rate, 7 bit data, 1 stop bit, and even parity for

Three modes of transmission are available: AUTO or POLL or MODBUS:

- In AUTO mode, the current PLC data string is transmitted at the end of analytical cycle without any request from the host computer.
- In POLL mode, the PP1 will transmit the latest PLC data string only when prompted by the host computer. This mode is recommended for use in a RS-485 network but can also be used for RS-232 connections.
- In MODBUS mode, the PP1 will transmit the latest PLC data string only when prompted by the host computer. See MODBUS uses manual for operational details.

Selecting the POLL mode on the MAIN SCREEN configures the PP1 to wait for a request before transmitting.

The master controller of the RS-485 network must first request data from the PP1 in order for receive the COM2 port output data.

The protocol to request data from the PP1 consists of:

- 1) A start of transmission character
- 2) Followed by the unit serial number
- 3) An end of transmission character.

or

STX character (02 hex)

SSSS

ETX character (03 hex)

The requesting command format follows: **<STX>SSSS<ETX>** , where SSSS = serial number of unit is used as the RS-485 ID.

6.3.2 COM 2 Port PLC Output Format

The COM 2 port output format consists of unit serial number, date and time stamp, error alarm, stream number, followed by impurity data. This data format is transmitted at the end of every run automatically if the COM 2 mode is set to Auto. This format uses a general protocol envelope. This envelope consists of a start of transmission character followed by data, and end of transmission character. The start of transmission indicator is a STX character (02 hex) and the end of transmission indicator is a ETX character (03 hex). The fields within each record are defined as follows:

<STX>SSS,YYYY-MM-DD,HH:MM:SS,E,X,N,AREA,CONC,N,AREA,CONC,.....N,AREA,CONC,<ETX>

- SSS = Serial Number of Unit. This is also used as the RS-485 unit ID
- YYYY-MM-DD = Date (2013-12-05)
- HH:MM:SS = Time (12:23:44 = 12 hours 23 minutes 44 seconds)
- E = General Error Alarm (0 = No Error , "1" = Error)
- X = Stream Number
- N = Name of Impurity (x6)
- AREA = Are of Impurity (x6)
- CONC = Concentration of Impurity (x6)

NOTE: The least significant decimal of the concentration resolution represents 0.1 ppb. (CONC reading of "1" represents a 0.1 ppb concentration level)

6.4 Modbus Protocol and Register Map for Peak Laboratories Devices

This section describes the Modbus protocol for Peak Laboratories PP1 analyzer. It is assumed that the reader is familiar with the Modbus protocol and serial communications in general.

The following rules define the protocol for information transfer between a Modbus MASTER device and the Modbus SLAVE. The MASTER initiates and controls all information transfer on the communications channel. A SLAVE device never initiates a communications sequence. All communications activity occurs in the form of "PACKETS." A packet is a serial string of 7-bit ASCII bytes. All PACKETS transmitted by the MASTER are REQUESTS. All PACKETS transmitted by a SLAVE device are RESPONSES. At most, one SLAVE can respond to a single request from a MASTER. The PP1 does not support broadcast request packets. The PP1 shall only support the ASCII mode Modbus protocol. The PP1 is only a SLAVE device. The PP1 does not conform to all stated requirements according to the "MODBUS Application Protocol Specification v1.1b3" and "MODBUS over Serial Line Specification and Implementation Guide V1.02.

6.4.1 Mode of Transmission

The PP1 only supports the ASCII mode Modbus protocol. The PP1 analyzer requires the serial communications channel to be configured to 9600 bps, 7 data bits, even parity, and one stop bit.

6.4.2 Modbus Packet Structure

Every Modbus packet consists of four fields:

- Slave Address Field
- Function Field
- Data Field
- Error Check Field (LRC Checksum)

Slave Address Field

The slave address field of a Modbus packet is one byte, two ASCII characters, in length and uniquely identifies the slave device involved in the transaction. Valid addresses will be in a range between 1 and 247. A slave device performs the command specified in the packet when it receives a request packet with the slave address field matching its own address. A response packet generated by the slave has the same value in the slave address field.

Function Field

The function field of a Modbus request packet is one byte in length and tells the addressed slave which function to perform. Similarly, the function field of a response packet tells the master what function the addressed slave has just performed. On page 5 lists the Modbus functions supported by the PP1 when acting as slave, please refer to "Table 2: Modbus Functions Supported by the PP1 as Slave".

Data Field

The data field of a Modbus request is variable length, and depends on the function. This field contains information required by the slave device to perform the command specified in a request packet, or data being passed back by the slave device in a response packet. Data in this field is contained in 16-bit registers. Registers are transmitted in the order of high-order byte first, low-order byte second.

Example:

- A 16-bit register contains the value 62BE Hex, four ASCII characters.
- High order byte = 62 Hex, two ASCII characters
- Low order byte = BE Hex two ASCII characters
- This register is transmitted in the order 6 2 B E, four ASCII characters.
- Error Check Field (LRC Checksum, two ASCII characters)

Error Check Field

In Modbus ASCII mode, an 8-bit Longitudinal Redundancy Check (LRC) algorithm is used to compute the checksum byte. The checksum field enables the receiving device to determine if a packet is corrupted with transmission errors. The sending device calculates the checksum, 8-bit value, based on every byte in the packet, using the LRC algorithm. The calculated value is inserted in the error check field. The receiving device performs the same calculation, without the error check field, on the entire packet it receives. The resulting value is compared to the error check field. Transmission errors are indicated when the calculated checksum does not equal the checksum stored in the incoming packet. The receiving device ignores a bad packet.

To calculate the LRC:

1. Add up all the data bytes in the message (before converting to ASCII and without the initial colon and final CR/LF).
2. Throw away any bits that carry over 8 bits.
3. Make the result negative (by twos complement) to get the LRC byte.
4. Convert from the 8-bit hex value to two ASCII characters.

6.4.3 Packet Communications:

The PP1 analyzer only supports two Modbus functions, Read Holding Registers and Preset Single Register. Function codes 03, 06 respectively.

Modbus Function	ModBus Description	Function Description
03	Read Holding Registers	Fetch the current value in one or more holding registers from the PP1.
06	Preset Single Register	Write a specific values into a single holding register to the PP1.

Table 1: Modbus Functions Supported by the PP1 as a Slave

Function 03: Read Holding Registers:

To read PP1 parameter values, a Master must send the PP1 a Read Holding Registers request packet. The Read Holding Registers request packet specifies a starting address and the number of 16-bit registers to read. The Modbus registers are numbered from 40001 to 49999. This Modbus numbering corresponds to an address map of 0 to 9998. The PP1 responds with a packet containing the values of the registers in the range defined in the request. See PP1 address map table to verify address range and content.

Read Holding Registers Packet Example:

The next example, the PP1's Modbus slave address is 92 and the master requests to read two parameters from the PP1. These two parameters are requested from Modbus registers 40013 and 40014. In accordance with the Modbus protocol, holding register 40013 is numbered as address 12 when requested.

Slave Address: 92 = 5C (Hex)

Modbus Function: 03 = 03 (Hex)

Starting Register: 12 = 000C (Hex)

Num of Register: 02 = 0002 (Hex)

LRC Checksum : -109 = 93 (Hex) = 39 33 (ASCII)

Request Packet: non-shaded background denotes the DATA field of the packet.

Slave Address		Modbus Function		Starting Register (40013)				Num of Registers (2)				LRC Checksum	
5	C	0	3	0	0	0	C	0	0	0	2	9	3

Response Packet:

Slave Address		Modbus Function		Byte Count		Register 1 Address (12)				Register2 Address (13)				LRC Checksum	
5	C	0	3	0	8	4	E	3	8	2	A	6	0	8	9

The Master device receives the above data packet in response to the request packet:
 Modbus Register 40013 = 4E38(Hex) = 20024
 Modbus Register 40014 = 2A60(Hex) = 10848

Function 06: Preset Single Register:

Preset a value into a single holding register. Place specific values into a single holding registers in the PP1. There is only one holding register that can be written to the PP1. See PP1 register map below.

Preset Single Register Packet Example:

In the next example, the PP1's Modbus slave address is 143 and the master requests to write the value of 3 to Modbus register (40004). In accordance with the Modbus protocol, register 40004 is numbered as address 3 when requested.

Slave Address: 143 = 8F (Hex)

Modbus Function: 06 = 06 (Hex)

Starting Register: 3 = 0003 (Hex)

Preset Data: 3 = 0002 (Hex)

Request Packet: non-shaded background denotes the DATA field of the packet.

Slave Address		Modbus Function		Starting Register (40004)				Preset Data(3)				LRC Checksum	
8	F	0	6	0	0	0	3	0	0	0	2	6	6

Response Packet:

Slave Address		Modbus Function		Starting Register (40004)				Preset Data(3)				LRC Checksum	
8	F	0	6	0	0	0	3	0	0	0	2	6	6

As you can see from the response packet above, the response is the echo of the requested packet.

6.4.4 Modbus Data Format:

The PP1 is capable of exporting all holding registers into the Modbus protocol.

The PP1 data available in a contiguous set of Modbus holding registers. The PP1 output registers are located in the Modbus register map from 40001 to 40099. The actual location is defined in table 3 below.

Six Data Formats Used in the PP1 Analyzer:

- UINT16 - 16 bit unsigned integer 0 to 65,535
- INT16 - 16 bit signed integer -32,768 to +32,767
- UINT32 - 32 bit unsigned integer 0 to 4,294,967,295
- INT32 - 32 bit signed integer -2,147,483,648 to +2,147,483,647
- Boolean - Packed Boolean Format
- 8CharString – 8 character string format

16-bit Integer Format:

Unsigned and signed 16 - bit integer formats are the simplest formats. Each PP1 address corresponds to one Modbus Holding Register. If the format is unsigned, the value range for the output registers is 0 to 65535. If the format is signed, the value range is -32767 to +32767.

32-bit Integer Format:

To accommodate values that can reach beyond the 16- bit range, the Modbus Slave provides 32 - bit integer format as an output option. In signed and unsigned 32 - bit integer formats, each PP1 address corresponds to two 16 - bit Modbus Holding Registers. A 32 - bit register represented in 32 - bit integer format is passed via communications as two 16 - bit registers:

High - Order Register = value / 65536

Low - Order Register = value modulus 65536

Example (Unsigned 32-bit):

Value 12345678 is passed in *unsigned* 32-bit integer format:

12345678 = 00BC614E Hex

High Order = 00BCHex

Low Order = 614EHex

Example (Signed 32-bit):

Value -12345678 is passed in *signed* 32-bit integer format:

12345678 = FF439EB2 Hex

High Order = FF43 Hex

Low Order = 9EB2 Hex

Packed Boolean Format:

Boolean registers are packed into a single Modbus register. The boolean outputs corresponds to one bit in the single output register of the module. The relationship is left to right: the first input register corresponds to the left - most bit in the 16 - bit output register, etc.

Example (Packed Boolean Format):

Six Boolean registers are linked to a single Modbus register, which is configured for Packed Boolean output format. If the first three are valued 'False', and the remaining three are valued 'True', the output register value is:

Register: 0001110000000000 Bin = 1C00 Hex

If the first input register became 'True', the output register value changes to:

Register: 1001110000000000 Bin = 9C00 Hex

8CharString Format:

The 8CharString or eight-character-string register format is four 16 - bit Modbus Holding Registers that hold a total of eight ASCII characters.

6.4.5 Broadcast Packets:

The purpose of a broadcast request packet is to allow all slave devices to receive the same request command from the master. A broadcast request packet is the same as a normal request packet, except the slave address field is set to zero (0). All Modbus slave devices receive and execute a broadcast request command, but no device will respond. The PP1 does not support broadcast request packets.

Exception Responses:

If a Modbus master device sends an invalid command to a PP1 or attempts to read an invalid holding register, an exception response is generated. The exception response follows the standard packet format. The high order bit of the function code in an exception response is set to 1. The data field of an exception response contains the exception error code. The table below describes the exception codes supported by the PP1 and the possible causes.

Return Function	Function Description
01	An invalid command is contained in the function field of the request packet. The PP1 only supports Modbus functions 03 and 06.
02	The address referenced in the data field is an invalid address for the specified function. This can also indicate that the registers requested are not within the valid register range of the PP1.
03	The value referenced in the data field is not allowed for the referenced register.
04	The LRC checksum is invalid.

Table 2: Exception Codes supported by the PP1

Exception Response Packet:

Slave Address		Modbus Function		Exception Code		LRC Checksum	
8	F	8	3	0	3	E	B

PP1 Modbus Registers Map:

The PP1 Modbus holding register map defines a set of parameters which are treated as holding registers, having addresses 4XXXX. According to the Modbus protocol, in response to a request for register 4XXXX of a particular slave device, the Modbus master reads register XXXX - 1 from the slave. For example, Modbus register 40054 corresponds to holding register 53. There are two classes of Modbus registers within the PP1, Modbus Output Registers, and External Control Registers. See Table 3 below.

PP1 Modbus Holding Registers Map

Modbus Register	PP1 Description	Format	Properties
40001	Unit serial number	UINT16	Read
40002	Main controller software version	UINT16	Read
40003	Detector type (RCP, FID, PDD, TCD)	UINT16	Read
40004	Run Mode (Idle, Single, Cycle, Stream)	UINT16	Read/Write
40005	Run Counter	UINT16	Read
40006	General Alarm Status	UINT16	Read
40007	FID Flame Out Alarm Status	UINT16	Read
40008	Heater Status Zone 1 (On/Off)	UINT16	Read
40009	Heater Set Point Zone 1	UINT16	Read
40010	Heater Temperature Zone 1	UINT16	Read
40011	Heater Status Zone 2 (On/Off)	UINT16	Read
40012	Heat Set Point Zone 2	UINT16	Read
40013	Heater Temperature Zone 2	UINT16	Read
40014-40039	Spare	-	-
40040	Compound #1 Name	8CharString	Read
40044	Compound #1 Area Counts	UINT32	Read
40046	Compound #1 Concentration (PPB/10)	UINT32	Read
40048-40049	Spare	-	-
40050	Compound #2 Name	8CharString	Read
40054	Compound #2 Area Counts	UINT32	Read
40056	Compound #2 Concentration (PPB/10)	UINT32	Read
40058-40059	Spare	-	-
40060	Compound #3 Name	8CharString	Read
40064	Compound #3 Area Counts	UINT32	Read
40066	Compound #3 Concentration (PPB/10)	UINT32	Read
40068-40069	Spare	-	-
40070	Compound #4 Name	8CharString	Read
40074	Compound #4 Area Counts	UINT32	Read
40076	Compound #4 Concentration (PPB/10)	UINT32	Read
40078-40079	Spare	-	-
40080	Compound #5 Name	8CharString	Read
40084	Compound #5 Area Counts	UINT32	Read
40086	Compound #5 Concentration (PPB/10)	UINT32	Read
40088-40089	Spare	-	-
40090	Compound #6 Name	8CharString	Read
40092	Compound #6 Area Counts	UINT32	Read
40096	Compound #6 Concentration (PPB/10)	UINT32	Read
40098-40099	Spare	-	-

Table 3: Modbus Holding Registers Supported by the PP1

6.5 Trend / Analog Outputs

There are seven independent analog outputs on the back of the PP1. The first six are tied to the six impurities as trend outputs and the last channel is a simple recorder output. These trend outputs are provided via 0 – 1.0 VDC screw terminals. Each impurity can be independently scaled by entering the desired range (e.g., 0 → 1000 = 0 to 1000 ppb = 0 – 1 VDC scale. At the end of each run, the analog trend signal is updated and holds until the end of the next run is completed. As an accessory, Peak can provide standard commercial VDC mA converters as needed.

6.6 Setting up a RS485 Network

The RS-485 specification supports two-wire half-duplex communications (only one unit may transmit at a time), but allows up to 32 users on a “party line” network.

NOTE: The RS-422 specification supports four-wire full-duplex communications (two units may transmit at a time) and also allows up to 32 users on a “party line” network.

The following information is provided to enable the users to install the Peak Performer 1 into a RS-485 network. This system diagram consists of the PP1, RS-232 to RS-485 converter, and master host controller.

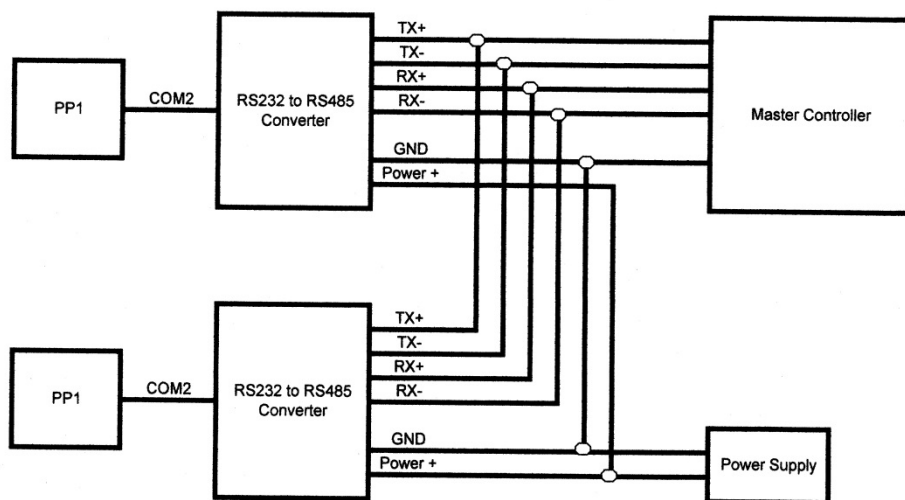


Figure 42. System configuration diagram with basic connections.

6.6.1 RS232 to RS485 Converters

There are numerous manufacturers of RS232 to RS485 converters and the selection of the converter must be determined by the system engineer / designer to meet unique requirements.

7.0 Chromatographic Principles

As an analytical system, the PP1 – TCD performs four primary functions:

- Sample injection
- Component separation
- Component analysis
- Integrated microprocessor system control with operator and data interface

7.1 Sample Injection

A standard air-actuated VICI injection valve is used to inject samples. Standard TCD analysis methods use stainless steel 1 cc sample loops.

An on-column injection adapter is available from Peak Labs if smaller sample volumes are required, using standard micro-volume or Luer type syringes.

7.2 Carrier Gas Purification

Chromatographic instrument detection limit is directly related to carrier gas purity. Improved carrier gas purity enables improved sensitivity.

Source	High Pressure Cylinder or Liquid Dewar
Gas Type	Ar or He
THC Concentration	< 10 ppm
CO, CO₂ Concentration	< 10 ppm
O₂ Concentration	< 10 ppm
H₂, H₂O Concentration	<10 ppm

7.3 Component Separation

Component separation is normally performed by isothermal packed analytical columns. Column materials vary by application and are noted in the test report documentation. Please contact Peak Laboratories for additional details specific to your application.

In general, applications samples are introduced into the instrument through the SAMPLE IN connection. A controlled portion of the sample gas is injected via the pneumatically actuated gas sampling valve into the carrier flow path.

The injected sample passes onto analytical column(s) inside the isothermal column oven (see figure 43) and moves over first (stripper) column, where the separation of the components begins. Analytes of interest are allowed to continue onto the second (analytical) column and the remaining balance of the sample is backflushed to vent. The analytes of interest are further separated on the second column and continue to the detector inlet.

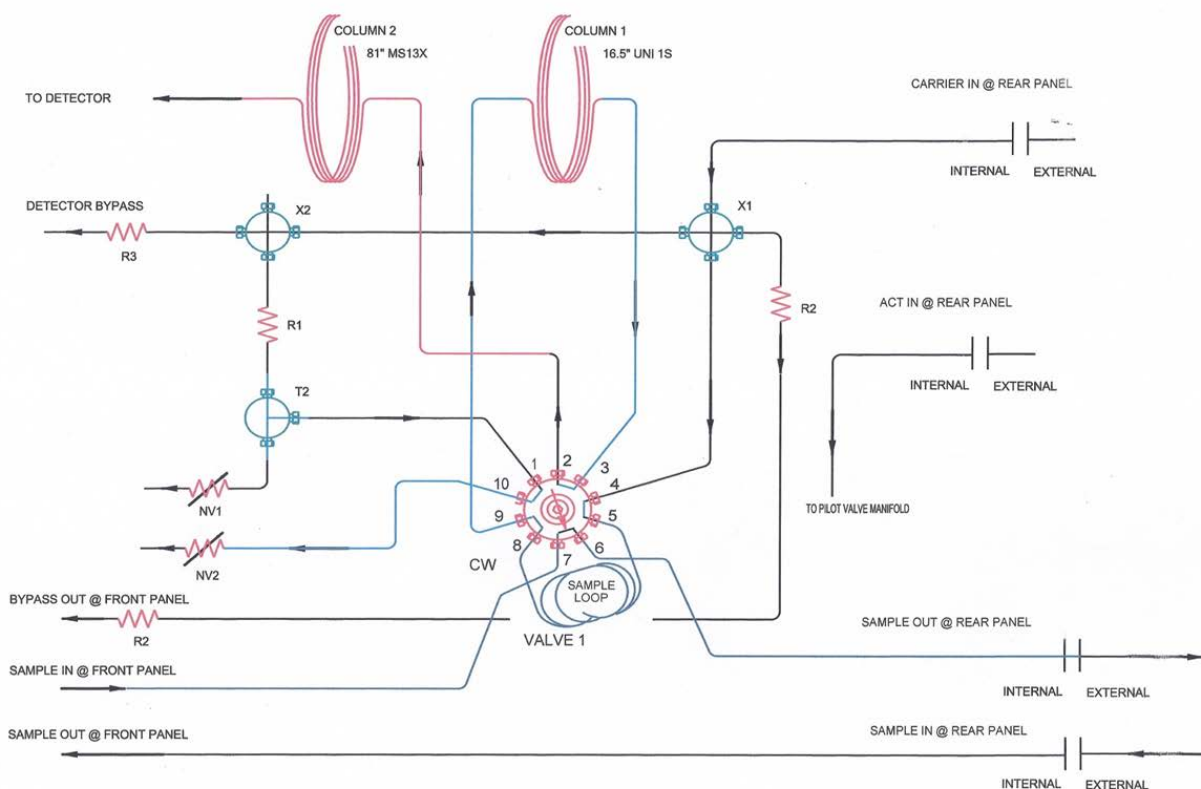


Figure 43. Typical TCD plumbing diagram

Immediately after injection, the entire contents of the sample loop are moved to the head of the column(s) and component separation begins.



Figure 44. Peak Development on the Analytical Column Soon After Injection

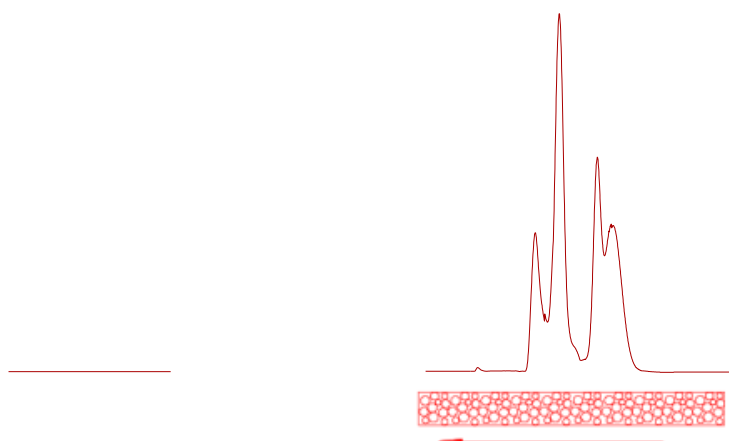


Figure 45. Peak Development Midway on the Analytical Column

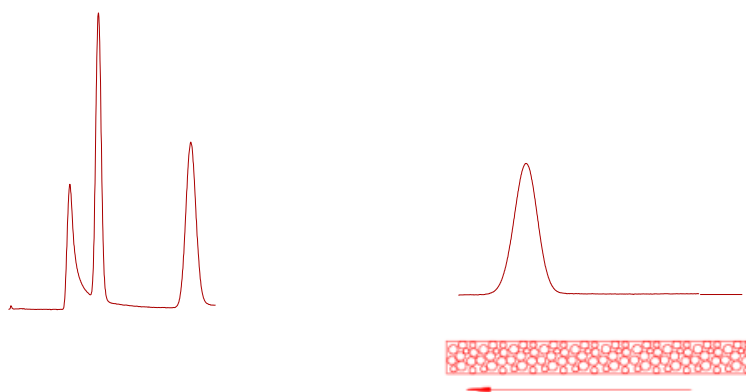


Figure 46. Peak Development Midway on the Analytical Column

Analytes of interest can be allowed to continue on to the detector as long as desired. In some cases, the resulting chromatogram can become quite long in duration. In other circumstances, it is not necessary to quantify all possible peaks, and the remaining balance of the sample is backflushed to vent or detector.

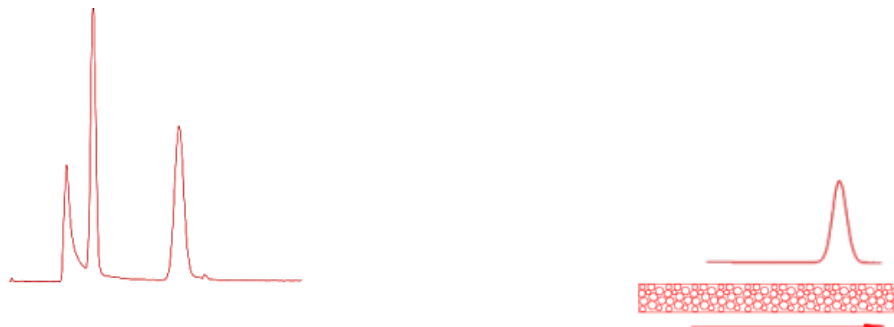


Figure 47. Peak Development on the Analytical Column After Backflush



Figure 48. The Completed Chromatogram

7.4 Flow Venting and Balancing

Several problems can arise in the development of the chromatogram:

- Presence of compounds that may be damaging to the detector
- Disproportionate peak sizes
- Length of analysis

Several models of Peak analyzers are equipped with one or two valve plumbing arrangements to handle these situations. Often a second valve is used to perform “cutting” actions at particular times during the chromatogram, redirecting the column effluent towards a safe flow controlled “vent” outlet.

During the “cutting” action, it is necessary to maintain pure gas flow to the detector. This pure gas flow is provided by an additional branch on the internal carrier gas main line feeding other ports on the valve.

Ultra-high purity flow control components can be bulky and expensive. To diminish installed cost of flow control components, Peak has chosen a different method based upon flow restrictors to create a flow divider. A typical flow divider is shown in Figure 49.

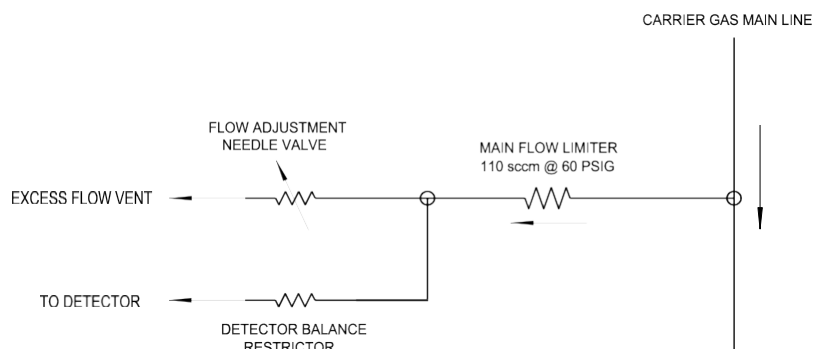


Figure 49. Flow Divider Schematic

The gas flow path to the detector is extremely clean and free of moving parts. Gas passing from the main flow limiting restrictor can be split to follow both gas paths. Studying the schematic, it can be seen that the detector flow can be adjusted by changing the amount of gas that is allowed to exit via the excess flow vent; as more gas exits the excess flow vent, less gas is available to flow towards the detector.

The controlling needle valves for detector flow rate and column venting are co-mounted in the chassis for convenience. The detector flow control needle valve is color-coded red.

Differences in DETECTOR flow during normal and venting valve positions can be expressed as shifts in detector baseline. Careful adjustment of the detector flow control needle valve will enable the operator to minimize baseline drift and create a smooth chromatogram.

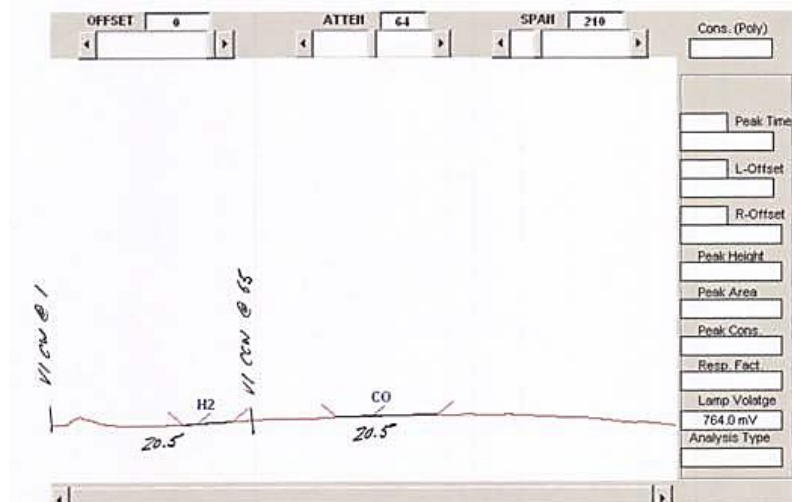


Figure 50. Properly Flow Balanced Baseline

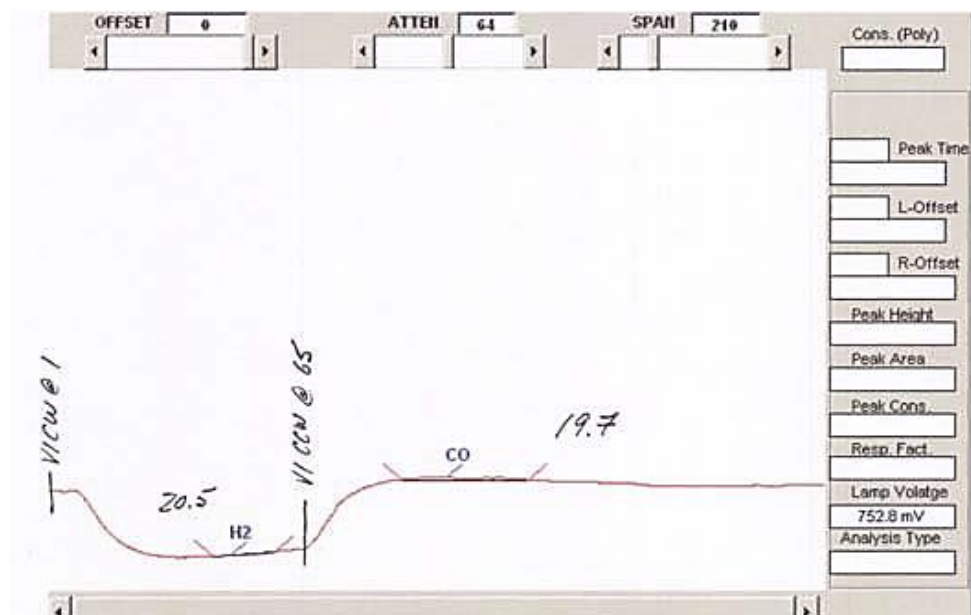


Figure 51. High TCD flow from NV1 (note rise after the V1 CCW command)

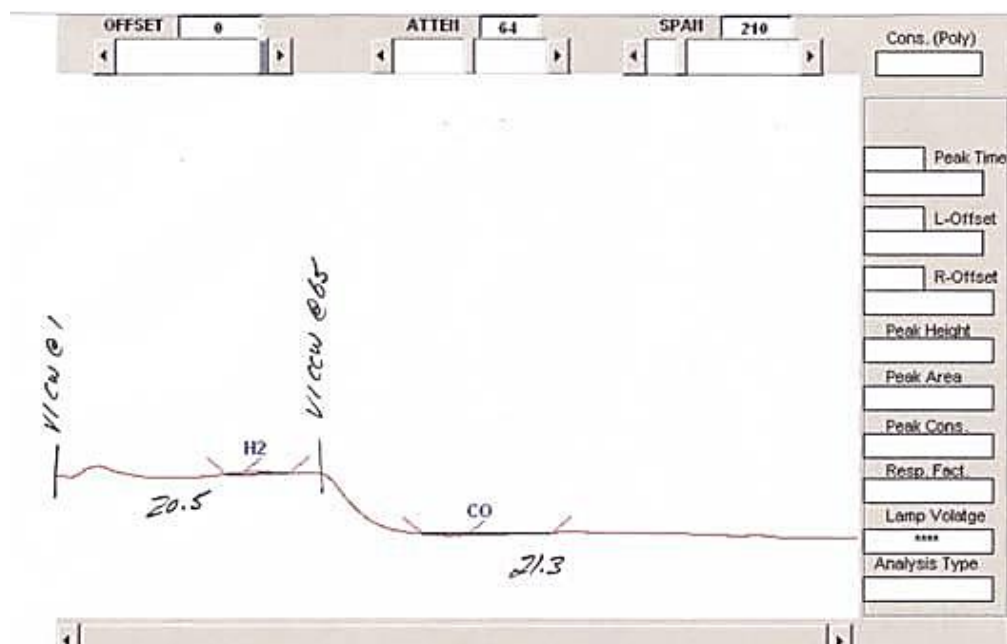


Figure 52. Low TCD flow from NV1 (note drop after the V1 CCW command)

Differences in COLUMN flow during normal and venting valve positions will be expressed as shifts in peak retention time compared to a chromatogram without any venting action.

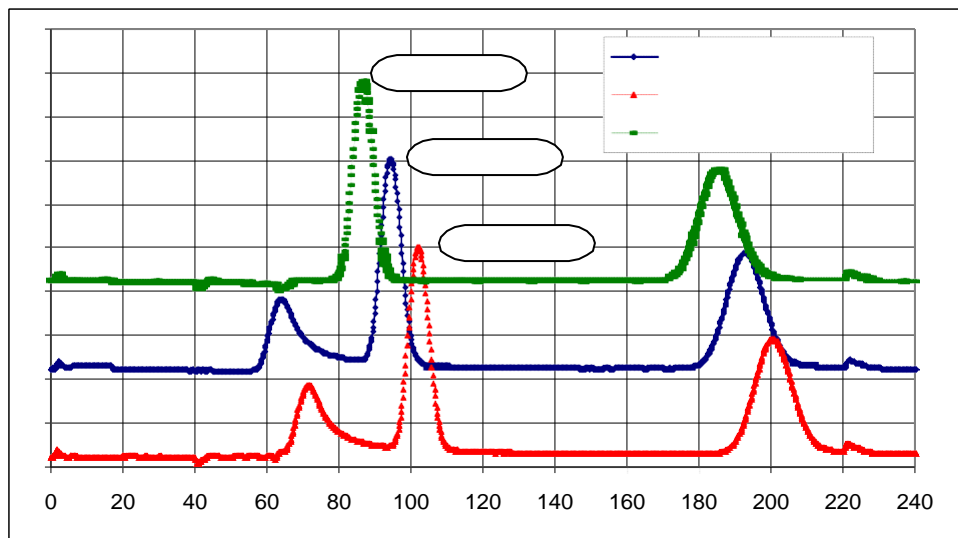


Figure 53. Peak Shifting as a Result of Improper Vent Needle Valve Flow

Comparison of retention times for un-vented versus vented chromatograms will reveal improper adjustment of the vent flow control needle valve.

Careful adjustment of the vent needle valve will enable the operator to minimize peak retention drift and create a consistent chromatogram.

7.5 Effects of Temperature and Flow on Analytical Columns

The normal chromatogram can be affected by changes in column flow and temperature. Be aware that changes of operating parameters from final test conditions can require parametric adjustments (i.e. peak detection parameters, response factors, valve timing) within the analyzer.

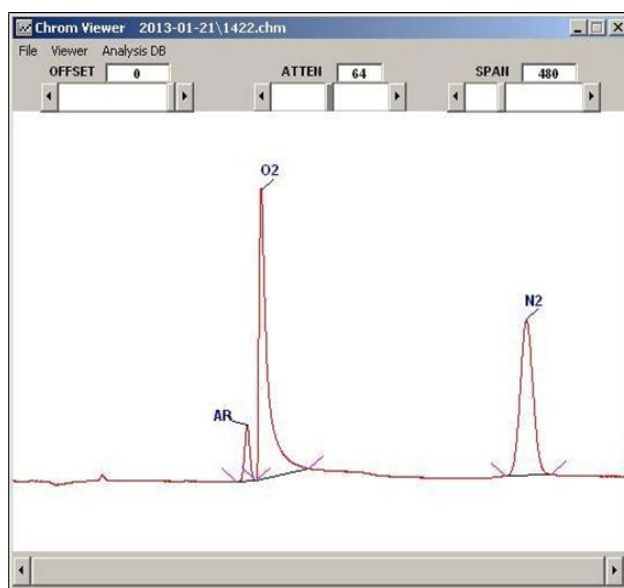


Figure 54. Standard chromatogram

The normal PP1 chromatogram (Figure 54) is affected by changes in column temperature as illustrated in comparison with the chromatogram of Figure 55.

In the circumstance of column temperature change, an increase of 30C has shifted the late eluting peak (carbon monoxide) forward. Early eluting peaks such as hydrogen are not significantly affected by changes in column temperature as interaction with column packing is minimal.

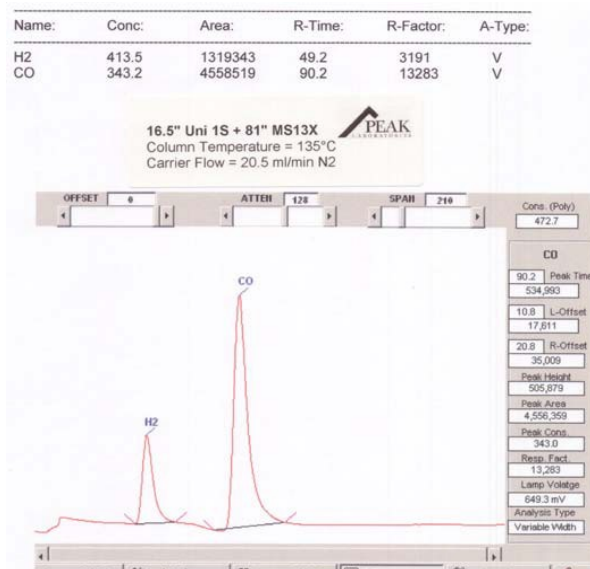


Figure 55. Effect of 30C Increase in Column Temperature

Column flow rate changes can directly affect peak retention times as peaks are pushed down the column faster or slower. Comparison of Figure 56 with Figure 54 shows the decreased peak retention times

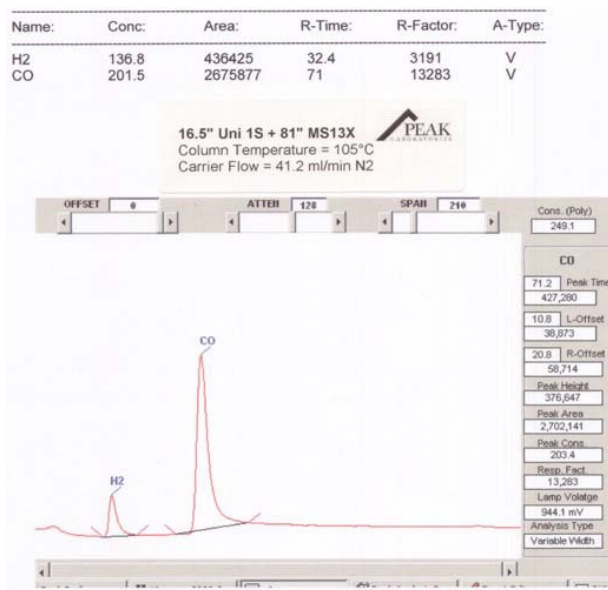


Figure 56. Effect of 2X Increase in Column Flow Rate

7.6 Peak Identification and Quantification

A typical chromatographic peak is generated by the detector's electronic signal as a function of time

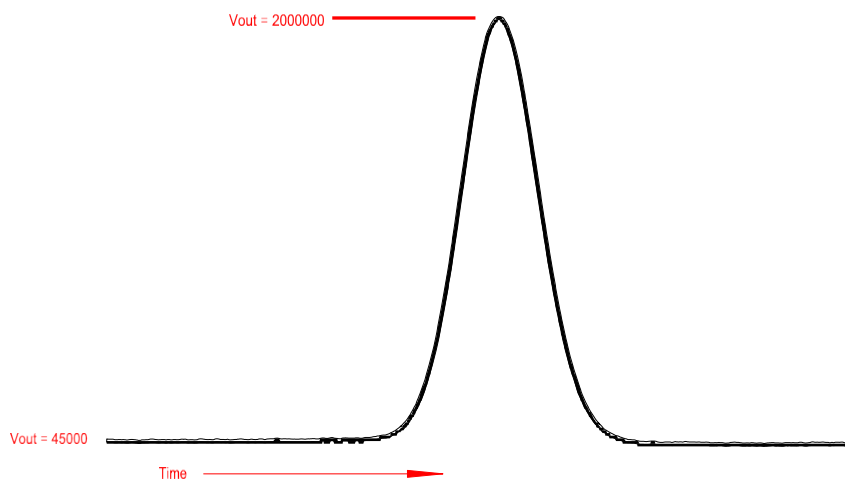


Figure 57. Typical Chromatographic Detector Signal Showing a Peak

Chromatographic peaks in the PP1 are measured by establishment of the “normal” stable detector signal, then determining the “added” electronic signal that was created by the peak passing through the detector.

Since the analyzer is monitoring the detector electronic signal over time, the “added” signal is measured in real units such as volt-sec. By virtue of a conversion factor (called the Response Factor), the volt-sec signal can be directly related to concentration units.

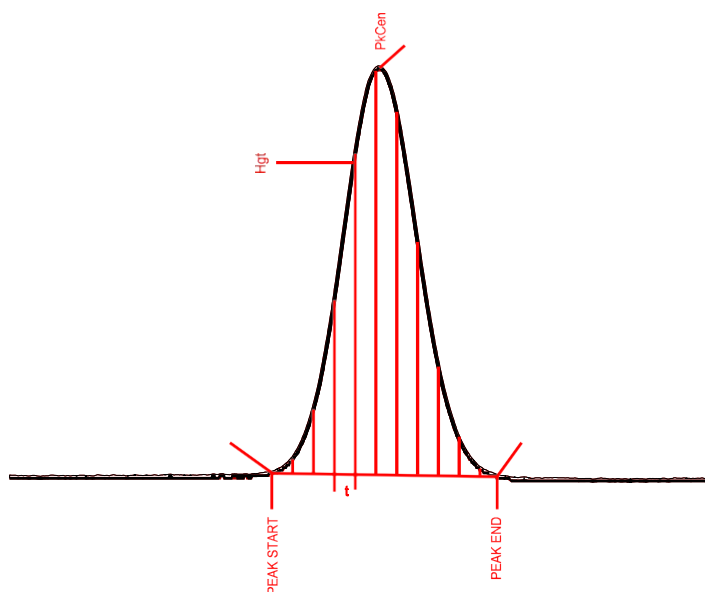


Figure 58. Illustration of Measurement of Chromatographic Peak Area

Note: The PP1 analyzer supports only area-based concentration determinations. Determinations based upon peak height are not currently available.

The first task in quantifying peaks is establishment of the “normal” stable detector signal (“baseline”). This is accomplished by determining the signal departure and return points for the peak event.

The parameters used when defining detector baseline signal and peaks are:

Analysis Database						
Name	PkCen	PkWin	PkHgt	LW	RW	Flt
*N2	226	10	3000	10	10	4

Analysis Averaging: 1

| ↑ | ↓ | Edit | Page

Figure 59. Typical PP1 Peak Analysis Database

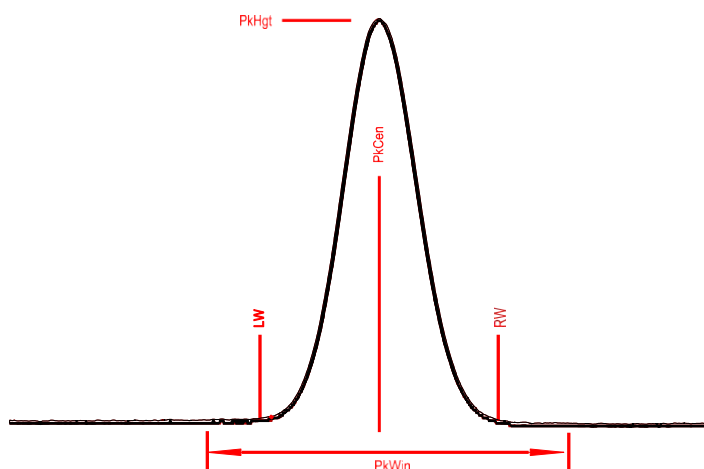


Figure 60. Illustration of Chromatographic Peak Parameters

Name = The peak identification tag, expressed as 5 alphanumeric characters (A-Z, 0-9).

PkCen = The normally expected peak retention time (in seconds), used for assigning the appropriate “Name”.

LW = The typical time span measured from the start of the peak’s baseline rise to the peak apex “PkCen”(in seconds). Minimum value is 3, maximum is 35.

RW = The typical time span measured from the peak’s apex “PkCen” to the end of the peak’s baseline decline (in seconds). Minimum value is 3, maximum is 35.

PkWin = Total tolerance window (in seconds) for assignment of a “Name” to a quantified chromatographic peak. The tolerance window is centered upon the “PkCen” value.
NOTE: The peak top detection algorithm processes only the data within this window.

PkHgt = A variable for establishing the cross-over point in peak detection modes. Chromatographic peaks higher than this parameter will be quantified using the “Variable” mode heights, peaks with height less than or equal to “PkHgt” will be quantified using the “ForceB” mode.

FIt = The convolution filter value expressing the overall peak shape. FIt = 2 is recommended for sharp, narrow peaks such as H₂, and FIt = 8 is recommended for broad peaks such as CO.

For good peak quantification, it is important to consistently determine when the peak begins (the positive increase in detector signal), when the signal maximum occurs, and when the peak ends (the return to stable detector signal).

PkCen, LW, and RW are used to establish starting and ending points for mathematical analysis of the baseline.

All modes of peak are calculation require locating the “peak top”, usually the point of maximum peak signal. In cases where the detector signal is strong, locating the peak maximum value is straightforward. Often the concentration regime of interest is the opposite case, where it is desirable to detect and extract the smallest possible peak signal.

Under these circumstances, it is advantageous to utilize signal processing techniques that magnify signal differences. Mathematical convolution is one such method. Convolution “rolls” two waveforms together as the product of data matrices, with the resulting new waveform having exaggerated characteristics.

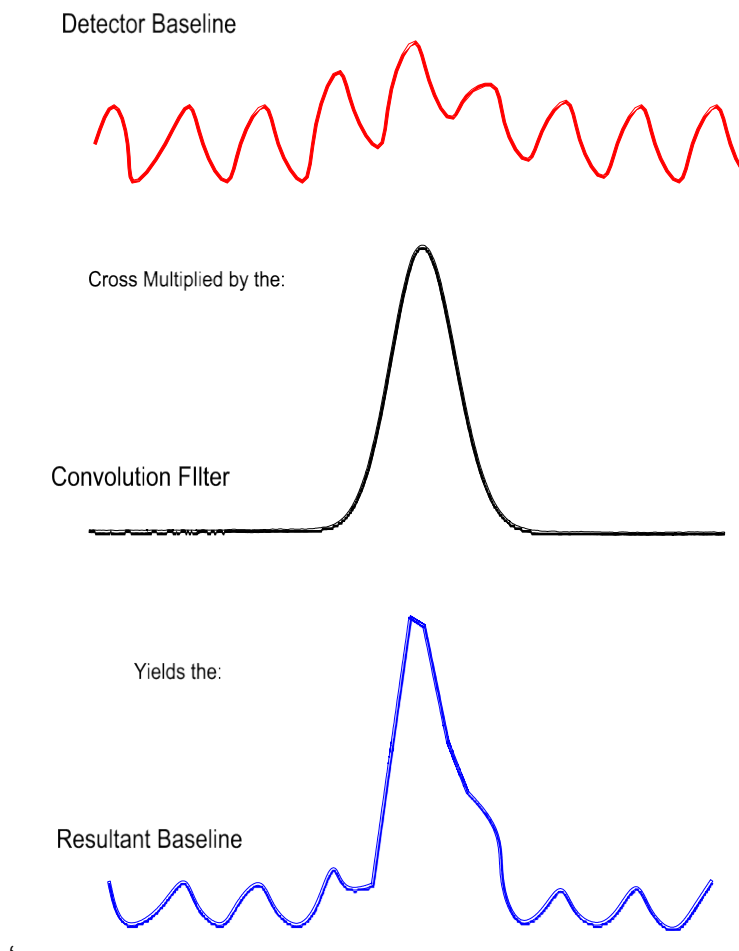


Figure 61. Convolution to Enhance Signal Characteristics

As seen in Figure 61, it is much simpler to identify the peak maximum on the convolved baseline.

Convolution functions work best when the convolving peak (filter) width closely matches the expected chromatographic peak width.

Chromatographic peak width can vary significantly depending upon column selection, detector characteristics, sample loop size, etc. Consequently, there are several different convolving

peak widths available in the PP1's mathematics package; Flt =1 is the narrowest filter width, Flt=8 is the widest filter width.

NOTE: Final determination of the optimum filter width is best performed by re-processing chromatograms representing typical peaks.

7.6.1 Peak Quantification: "Baseline" Mode

When the PkHgt variable is set to zero, the Baseline Mode of detector baseline is enabled. In this mode, the normal detector baseline is calculated very simply:

Peak start time = [PkCen – LW])

Peak end time = [PkCen + RW])

All other possible aspects of peak start and end time determinations are disabled and the normal detector baseline is rigidly established only by the Peak Database parameters.

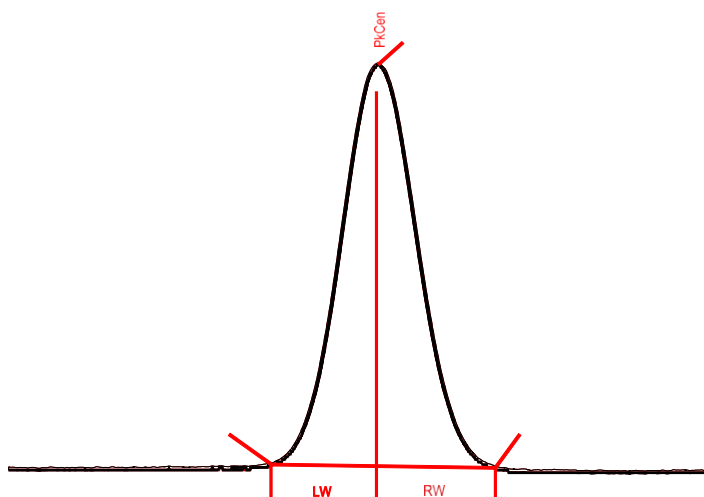


Figure 62. Baseline as Established in Baseline Mode

Peaks quantified using this method are denoted with an "B" in the last column of the data report

7.6.2 Peak Quantification: "Fixed Width" Mode

The Fixed Width Mode is enabled when the maximum detector peak signal value is less than the PkHgt parameter value, but is greater than zero.

The time associated with the peak maximum signal is always calculated first for Forced Baseline and Variable modes by using the convolution filter to find the maximum detector signal point.

The time associated with the peak maximum signal must be within the time regime of interest (= [PkCen +/- ½ PkWin]) for any peak quantification to occur.

In Fixed Width Mode, the normal detector baseline is calculated based upon the peak maximum signal time value, rather than the PkCen value.

Peak start time = [Peak signal maximum time value – LW]

Peak end time = [Peak signal maximum time value + RW]

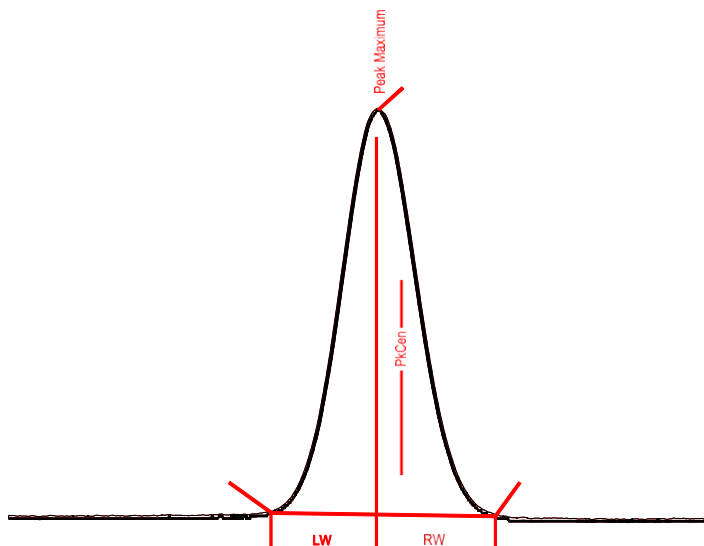


Figure 63. Baseline as Established in Fixed Width Mode

Peaks are quantified using this method are denoted with a “F” in the last column of the data report.

7.6.3 Peak Quantification: “Variable” Mode

If the peak maximum signal value exceeds the PkHgt value, the Variable Mode is enabled for determination of peak start and end points.

The time associated with the peak maximum signal is always calculated first for Forced Baseline and Variable modes by using the convolution filter to find the maximum detector signal point.

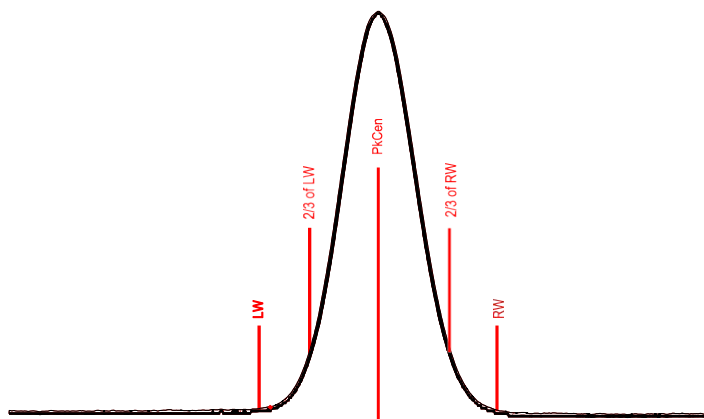


Figure 64. Initial parameters for Baseline Determination

Initially, the mathematics package makes a trial evaluation of peak area using a fixed percentage of the LW and RW values as initial baseline points, illustrated as the hypothetical baseline between points “2/3 of LW” and “2/3 of RW” in Figure 64. The hypothetical calculation of the peak start time is labeled “LW” in Figure 65.

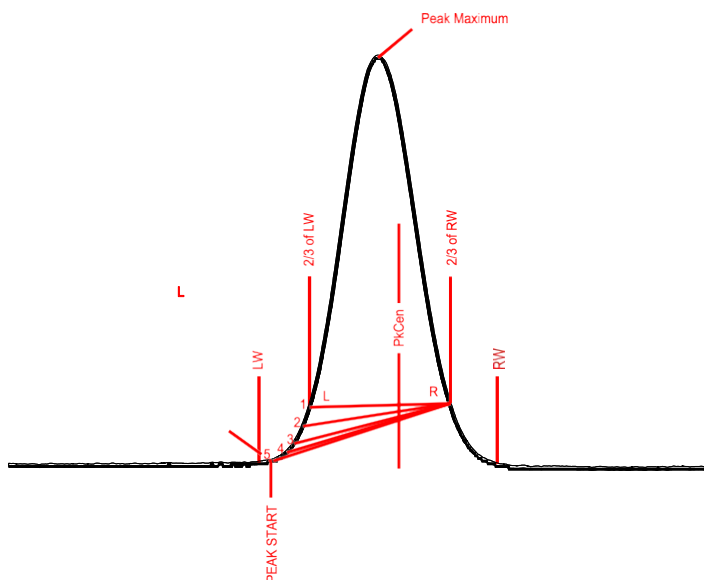


Figure 65. Progression of Peak Start Time Calculation

After calculating the hypothetical peak area based upon the L and R values, the mathematics package makes another peak area calculation, using a slightly earlier time value. The peak areas of the two calculations are compared, and if the peak area increases by more than +0.5%, the mathematics package repeats the process. The net effect is illustrated in Figure 65, where multiple hypothetical peak start points were evaluated. The final value is represented by point "5", when the change in peak area became less than +0.5%.

The trial evaluation of peak area continues for the peak ending point, again using a fixed percentage of the LW and RW values as initial baseline points. This is illustrated as the hypothetical baseline between points “L” and “R” in Figure 66. The first hypothetical calculation of the peak end time is labeled “1” in Figure 66.

After calculating the hypothetical peak area, the mathematics package makes another peak area calculation, using a slightly later time value. The peak areas of the two calculations are compared, if the peak area increases by more than +0.5%, the mathematics package repeats the process. The net effect is illustrated in Figure 66, where multiple hypothetical peak start points were evaluated. The final value is represented by point “6”, when the change in peak area became less than +0.5%.

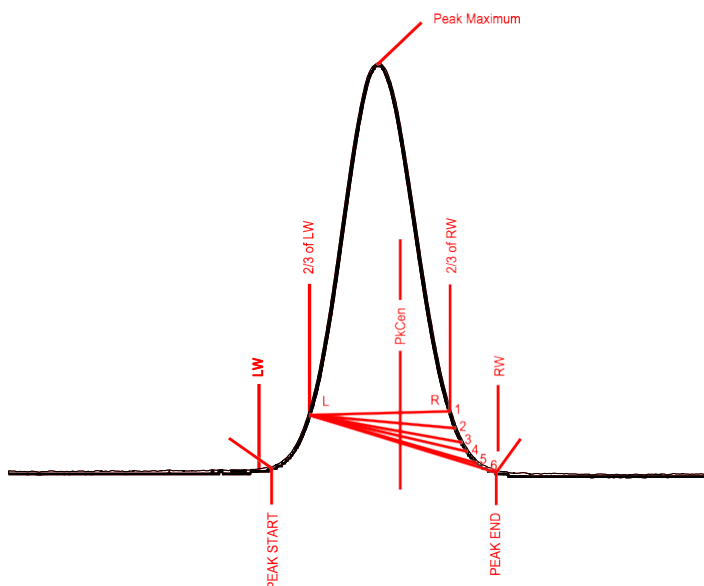


Figure 66. Progression of Peak End Time Calculation

Using the newly established peak start and end times, the mathematics package determines the “normal” detector baseline, as illustrated in Figure 67.

The signal contribution due to the peak passing through the detector can then be quantified by area calculation using time slices and signal differentials above the newly established baseline.

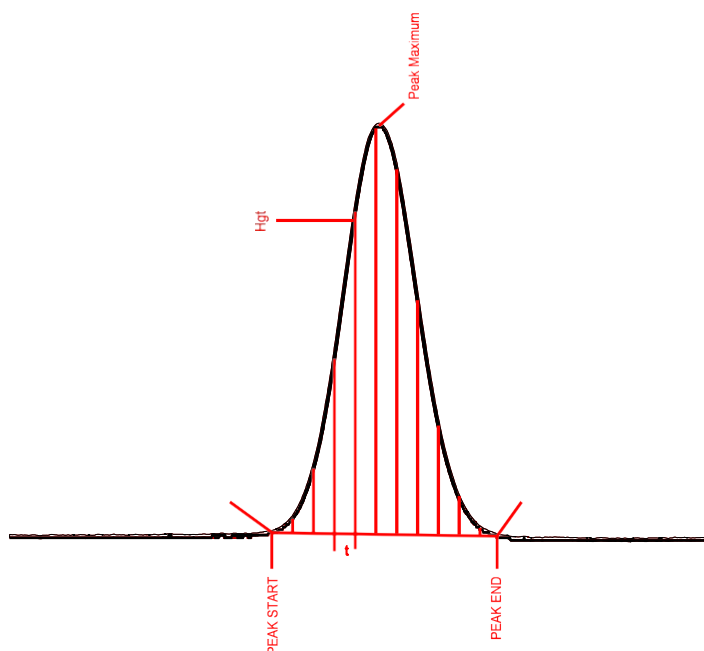


Figure 67. Peak Area Calculation

The formula for peak area calculation can be expressed as:

$$\text{Peak Area} = (\frac{1}{2} \times [\text{Hgt}_1 + \text{Hgt}_2] \times t_1) + (\frac{1}{2} \times [\text{Hgt}_2 + \text{Hgt}_3] \times t_2) + (\frac{1}{2} \times [\text{Hgt}_3 + \text{Hgt}_4] \times t_3) + \dots$$

where:

Hgt_x = the signal differential between the instantaneous baseline level and the peak signal value

t_x = the time difference between data points

Using a constant value for **t_x**, the peak area calculation simplifies to:

$$\text{Peak Area} = t_x \times [\text{Hgt}_1 + \text{Hgt}_2 + \text{Hgt}_3 + \text{Hgt}_4 + \dots]$$

Peaks quantified using the Variable method is denoted with a "V" in the last column of the data report.

8.0 Thermal Conductivity Detector (TCD)



**** CAUTION **** Do not begin detector heating without carrier gas flow. Damage to the detector may result.



**** CAUTION **** Never energize heaters of the instrument unless carrier gas flowing. Be sure to de-energize the filament before conditioning columns.



**** CAUTION **** The maximum detector operating temperatures for the TCD is 200°C.

The TCD detector does not normally require conditioning treatment, however it does require 20 to 24 hours at normal operating temperature to achieve optimum stable performance.

Normal detector operating temperature for the TCD is 60°C, although operating temperatures may vary with the specific application.

The table below illustrates some typical compounds and detection capabilities of the TCD.

Helium Carrier

Detected Compound	Typical Detection Limit
Hydrogen	10 ppm
Hydrocarbons	10 ppm
Argon	10 ppm
Nitrogen	10 ppm
Oxygen	10 ppm
Other Inerts	10 ppm

Figure 68. TCD Detection Capabilities

8.1 Principle of Operation of Thermal Conductivity Detection

Thermal conductivity detection relies on two fundamental characteristics of materials and gases:

- Transfer of heat (expressed as the coefficient of thermal conductivity)
- Changes in electrical resistance with temperature (expressed as the temperature coefficient of resistance).

In its simplest form, a thermal conductivity detector senses the heat loss from a heated electrical filament as a flowing gas passes over its surface. A usable detection signal can be generated since different gaseous compounds will conduct heat away from the filament at differing rates based upon the coefficient of thermal conductivity.

Gaseous Compound	Coefficient of Thermal Conductivity (<i>k</i>) W/m-°C
Hydrogen	0.168
Helium	0.124
Neon	0.046
Methane	0.030
Nitrogen / Oxygen / Air	0.024
Argon	0.016
Carbon Dioxide	0.0146
Xenon	0.0051

Figure 69. Coefficient of Thermal Conductivity

The heat loss from the filament can be measured by indirect methods such as a downstream temperature sensor, or by direct methods using the heated filament itself.

Peak Laboratories' TCD uses a direct measurement method, taking advantage of the natural variability of the filament resistance with temperature change.

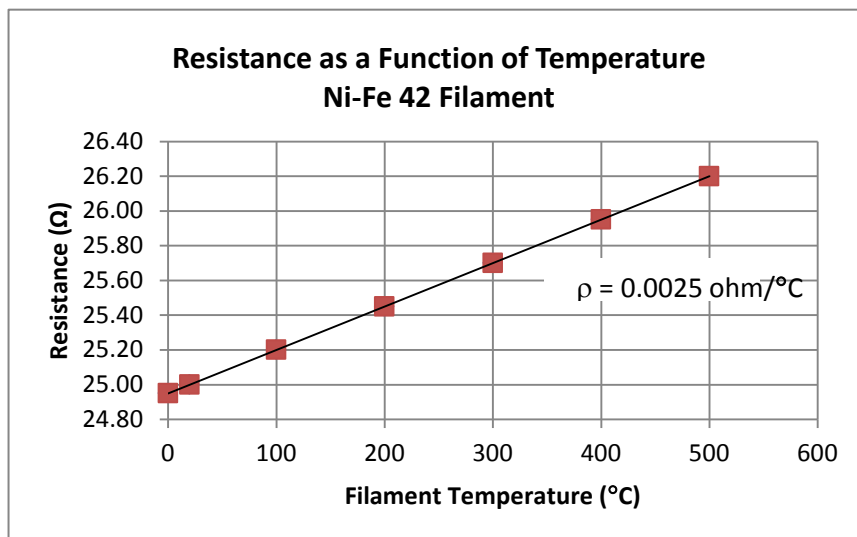


Figure 70. Resistance as a Function of Temperature

Temperature stability is critical for TCD operation, so the filament of the TCD is mounted in a precision controlled heated environment. When carrier gas flow is established through the TCD and the detector is energized, the filament itself is gently heated to a modestly higher temperature with an extremely stable constant current source. The applied voltage required to maintain the preset filament electrical current is monitored and amplified, providing the baseline detector signal.

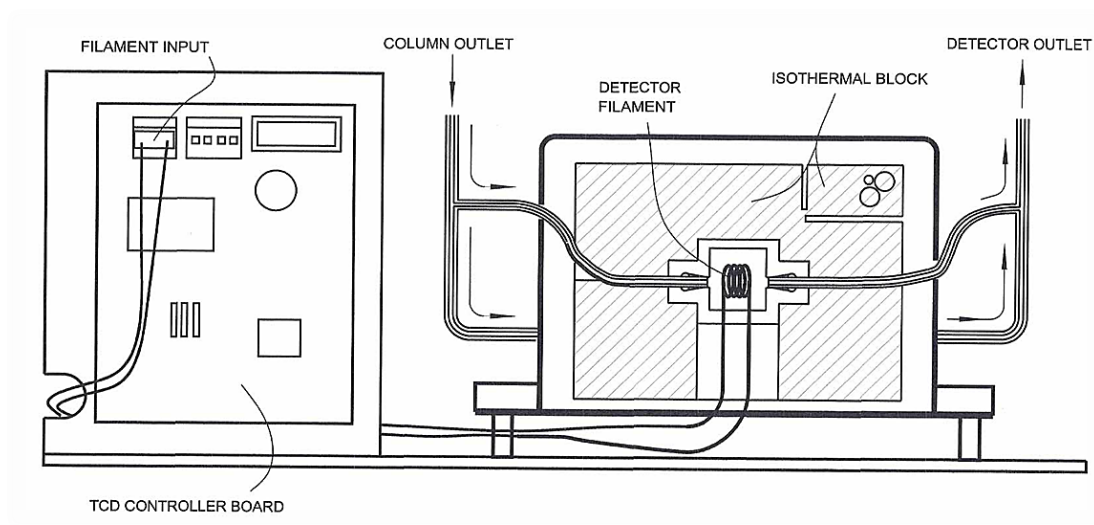


Figure 71. Standard TCD Detector

Once steady-state heat transfer from the filament to the flowing carrier gas stream is established, the filament resistance and the applied voltage stabilize.

In accordance with Ohm's Law, $V_{\text{(volts)}} = I_{\text{(amperes)}} \times R_{\text{(ohms)}}$, if the current value is held constant, changes in the resistivity of the filament will be reflected as changes in applied voltage.

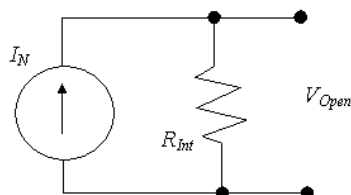


Figure 72. Simple constant current source

If a gaseous compound with a thermal conductivity different than the carrier gas flows past the filament, the steady-state heat transfer from the filament will be disturbed. Gaseous compounds with greater thermal conductivity than the carrier gas will extract more heat from the filament's surface creating a cooling effect. Conversely, gaseous compounds with lesser thermal conductivity will create an insulating effect at the surface of the filament, raising filament temperature.

As seen from Figures 70 and 71, these will create a small change in filament resistance, which in turn creates a small change in the applied voltage. The applied voltage is then amplified and digitized to create a chromatogram.

9.0 Analyzer, Maintenance and Service Procedures

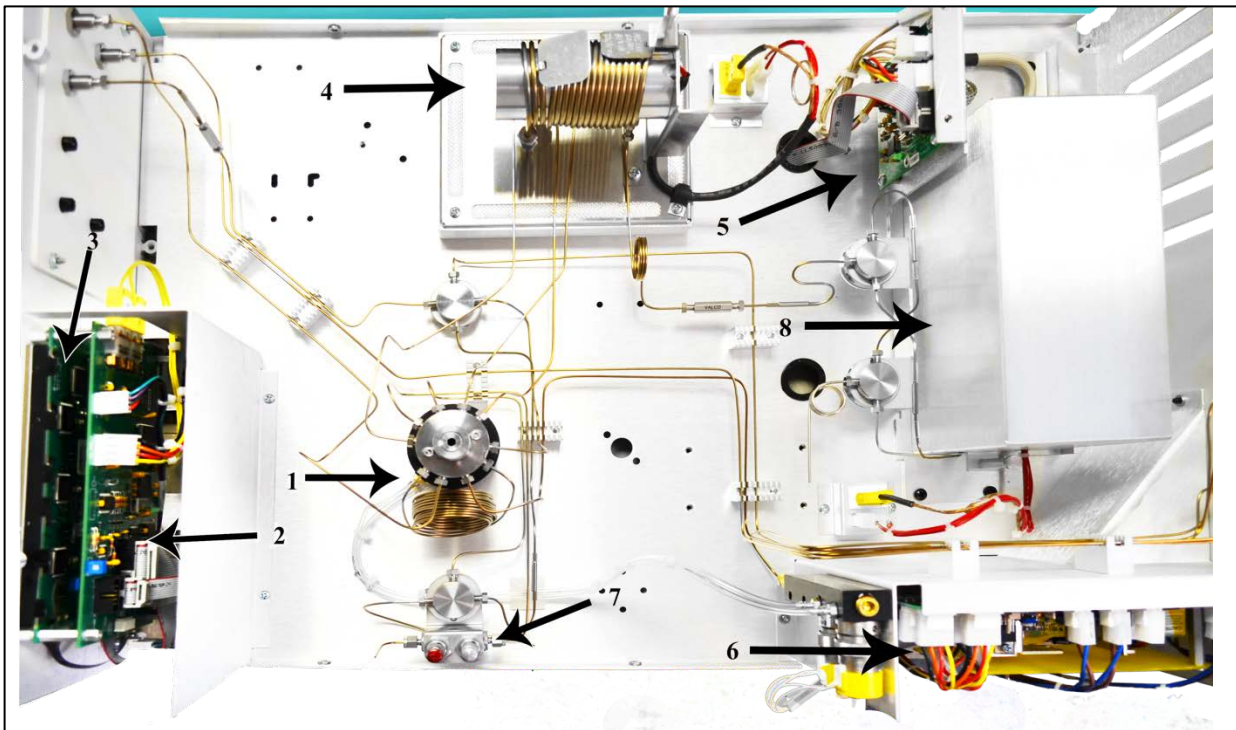


Figure 73. Chassis Overview

PEAK PERFORMER 1 TCD COMPONENT LAYOUT

1. VALVE 1 (10 Port)
2. MCPU PCB
3. LCD / TOUCHSCREEN ASSEMBLY
4. COLUMN OVEN without COVER
5. TCD CONTROLLER CARD
6. REAR PANEL CARD (INCLUDES MAIN DC POWER SUPPLY)
7. FLOW CONTROL NEEDLE VALVE ASSEMBLY
8. TCD DETECTOR WITH COVER

9.1 Peak Performer 1 System Block Diagram

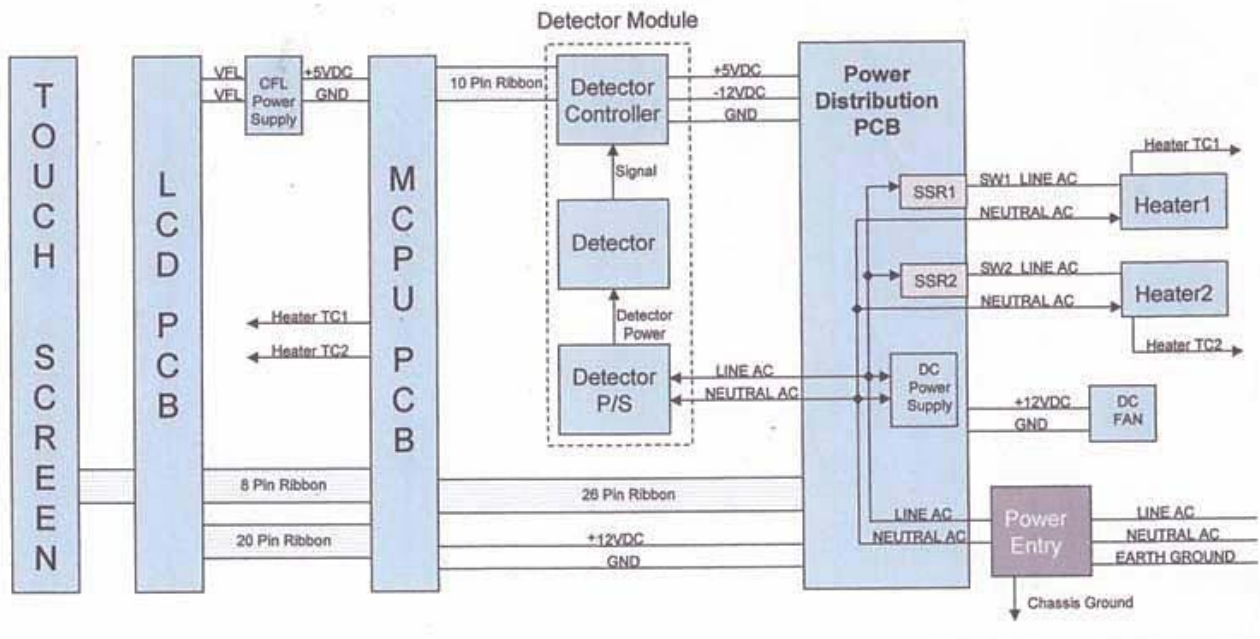


Figure 74. PP1 Functional Block and Interconnect Diagram

9.2 Monitor Regulator Pressures

Carrier supply settings should not vary. Carrier and Actuator Air settings may shift when supply cylinders run low due to “regulator creep” and should be checked periodically to maintain consistent gas flows and peak retention times.

Varying carrier pressures can cause peaks to shift out of set retention times and not be flagged properly.

Critical operating parameters are listed on the top page of the Final Test Data shipped with the analyzer. Peak Labs recommends the user create a logbook of critical parameters on a routine basis as an aid to troubleshooting.

Suggested parameters for frequent logging are:

- Analyzer S/N:
- Date and Time:
- Carrier Gas Type:
- Carrier Gas Supply Pressure:
- Carrier Gas Flow at NV1 Outlet:
- Carrier Gas Flow at NV2 Outlet, if applicable:
- Carrier Gas Flow at Front Panel Bypass Out Port:
- Carrier Gas Flow at NV1 Outlet (V1 CCW), if applicable:
- Carrier Gas Flow at NV2 Outlet: (V1 CCW), if applicable
- Actuator Gas Supply Pressure:
- Sample Gas Type:
- Sample Gas Flow at Rear Panel Sample Out Port:
- Column Setpoint and Measured Temperature:
- Detector Setpoint and Measured Temperature:
- TCD Vout Signal (after zeroing):
- Rfactor for compounds 1-6 (from Calibration Database screen)
- PkCen for compounds 1-6 (from Analysis Database screen)
- Quality Number

Other critical parameters that should be verified occasionally:

- Actuator Gas Type:
- Port 1 Protocol Setting
- Port 2 Protocol Setting
- Calibration Database Screen Parameters (all)
- Analog Interface Screen Parameters (all)
- Analysis Database Screen Parameters (all)
- Event Program Editor Screen Parameters (for all active Event Programs)

9.3 Verifying Valve Integrity (Carrier Blank, Zero Gas)

In addition to periodic calibration checks, the integrity of the rotary valve seal must be tested.

Determination of valve seal quality is performed by analyzing known purified gas as a sample.

- 1) Connect the front panel SAMPLE IN port to the front panel BYPASS OUT port. This will begin flow of purified gas through the sample loop.
- 2) Wait 10 minutes to purge the analyzer completely of any residual gas samples.
- 3) Initiate a standard sampling run in CYCLE mode.
- 4) Collect 5 sample runs with Peak Viewer.
- 5) Average the concentration data of the 5 runs for each peak of interest.
- 6) Average readings higher than 1 ppb for any compound indicates a loss of valve seal integrity and possibly suggests valve head replacement.

9.4 Adjusting the LCD Contrast

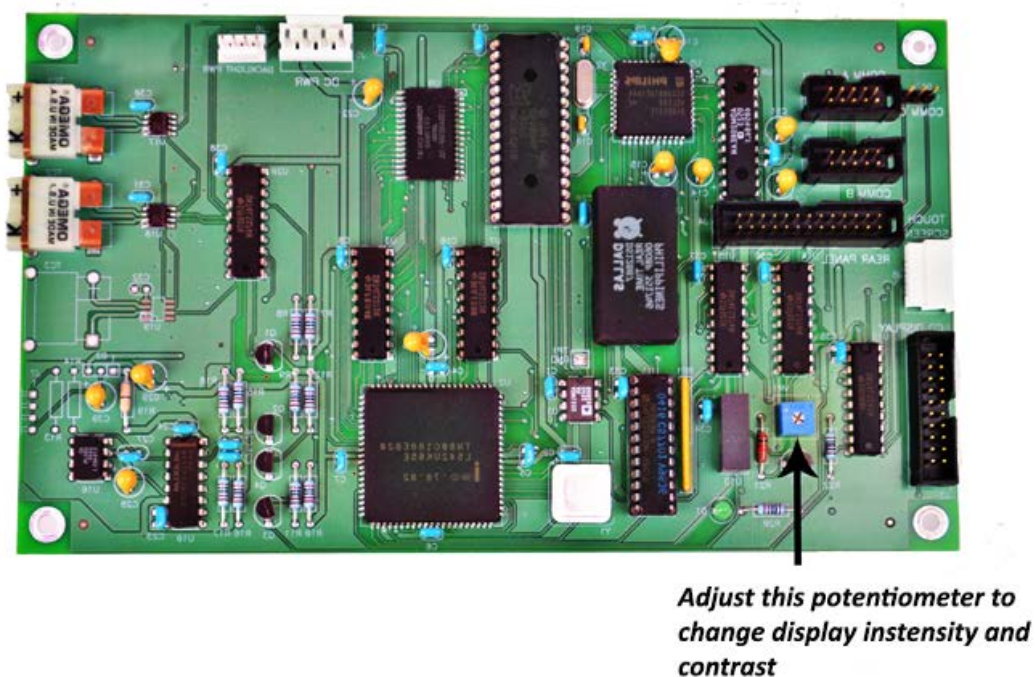


Figure 75a. View of the MCU Showing the LCD Contrast Potentiometer
(MCPU adjustments for software versions 1-3)

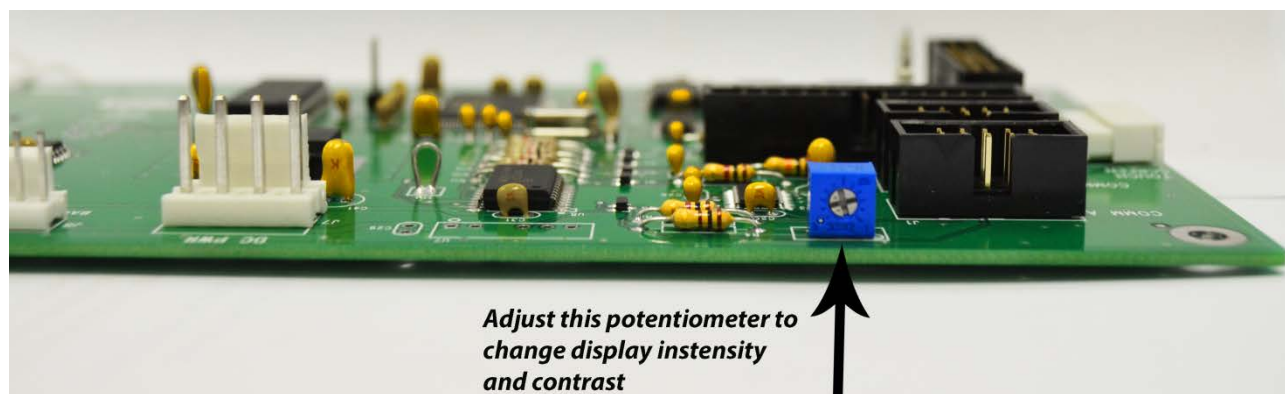


Figure 75b. View of the MCU Showing the LCD Contrast Potentiometer

9.5 ALTERING THE ANALYTICAL RANGE OF THE INSTRUMENT

Sample loops are very easy to customize.

- 1) Determine the desired loop volume. Analytical range is directly proportional to sample loop size
- 2) Select tubing (tube internal diameter) to use. In general, use the largest ID tube possible to avoid flow restrictions.

Guidelines:

0.040" ID for 3/8 cc to 2 cc loops
 0.030" ID for 3/16 cc to 1/2 cc loops
 0.020" ID for 50 uL to 1/4 cc loops

Note: the shortest practical length of tubing for a sample loop is about 8 inches.

Note: Do not use tubing less than 0.016" ID as it acts like a flow restrictor, instead of a sample loop

- 3) The correct length for the loop tubing is calculated as follows:

Internal volume (V in cc's) = [3.142 x ID ^ 2 (in inches)] x L (in inches) x 16.38 / 4 re-
 arranging:

$$L \text{ (in inches)} = 0.777 \times [V \text{ (in cc's)}] / [ID ^ 2 \text{ (in inches)}]$$

an example for a 1 cc loop:

$$L \text{ (inches)} = 0.777 \times 1 \text{ cc} / (0.040" ^ 2) = 49 \text{ inches tube length}$$

- 4) Measure to length, cut, and install VICI fittings
- 5) Blow the tube out with compressed air to remove any particles

9.6 Testing the Heater Resistance



Figure 76. Heater Resistance for 90 – 120 VAC AC Service

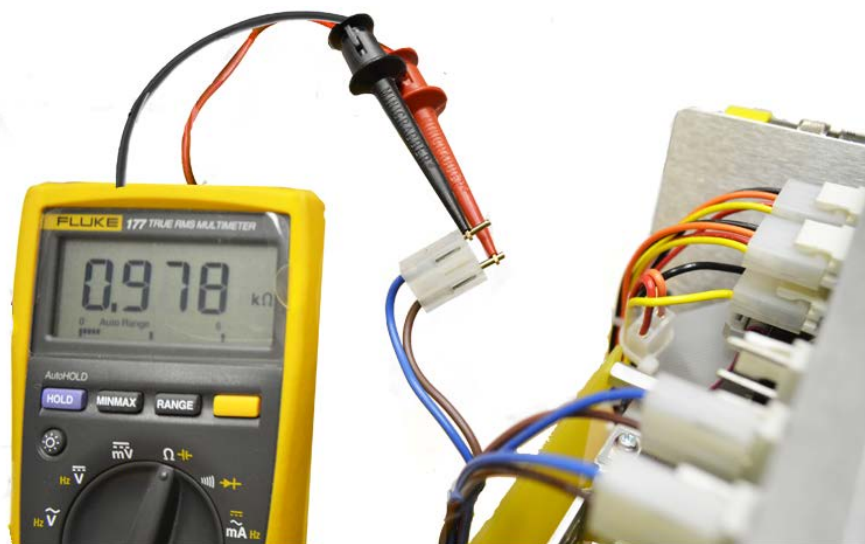


Figure 77. Heater Resistance for 200-240 VAC AC Service

9.7 Main DC Power Supply Replacement

Purpose:

Replace / upgrade Pihong PSA4531 DC Supply with Power-One BLP55-3000 DC Power Supply

Procedure:

Power off analyzer – if gas umbilical tubing is used, it is not necessary to cool down the analyzer.

Maintain normal gas supplies.

Disconnect AC power cord.

Remove top cover.

Locate the main DC power supply in the right rear portion of the chassis –see Figure 73.

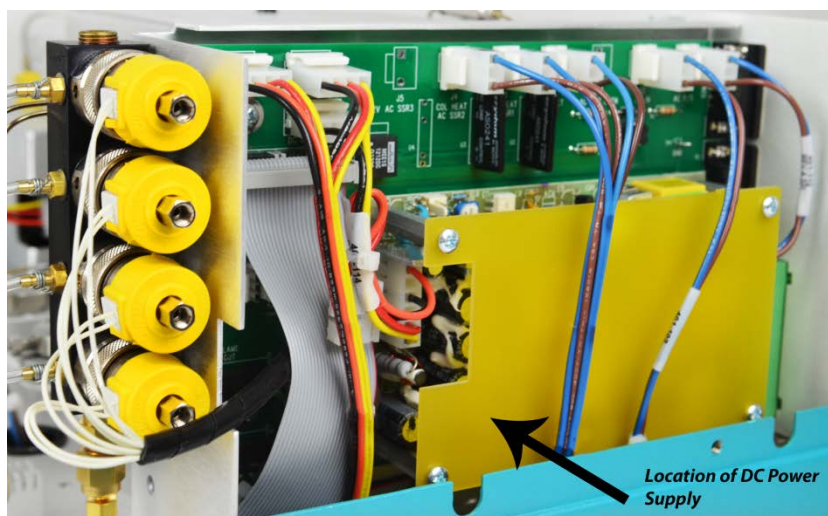


Figure 78a – Main DC Power Supply Location

Remove the 4 each 6-32 screws that hold down the plastic protective shield to expose the DC power supply – see Figure 75a.

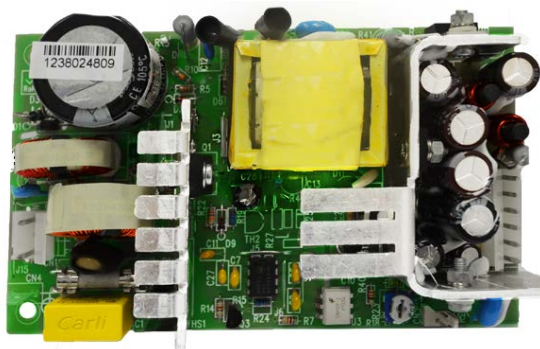


Figure 78b – Main DC Power



Figure 79. Main DC Power Supply with Shield Removed

Unscrew the 4 each 1/4" standoffs that retain the DC power supply to the rear panel circuit board. Disconnect only the two cable assemblies that attach the DC power supply to the rear panel circuit board – it is recommended to disconnect the cables at the circuit board connector end. Remove the DC power supply by sliding it under the overlaid cable assemblies – do not remove the other cable assemblies – see Figure 80.

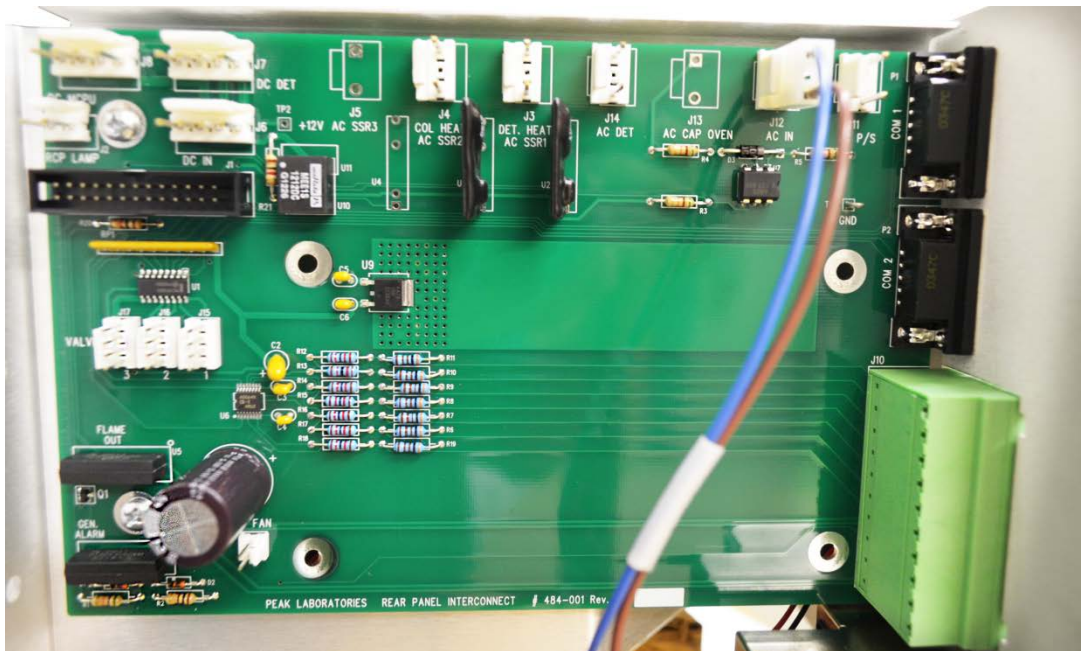


Figure 80. Main DC Power Supply Fully Removed

Lay the old DC power supply on the workbench aside the new Power-one DC power supply. Orient such that the 6 pin and 2 pin connectors are aligned – see Figure 79.


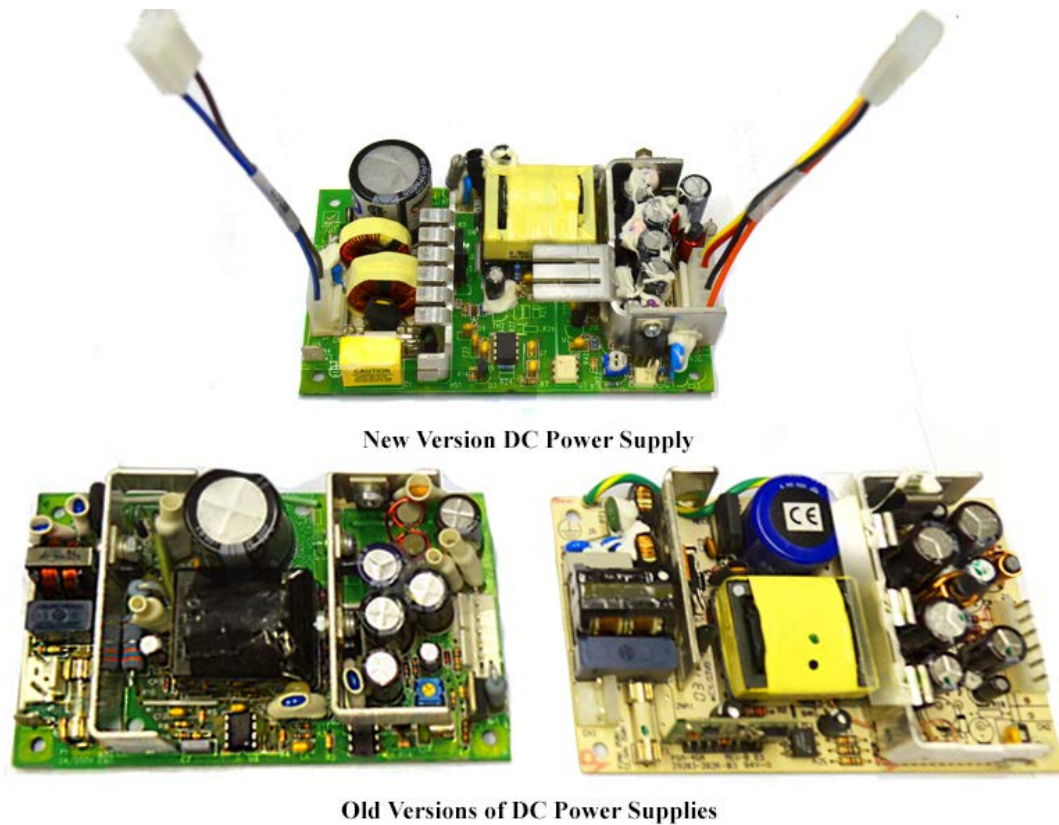


Figure 81. Correct Alignment and Cable Transfer from Old to New DC Power Supply

Transfer the two cable assemblies from the old power supply to the new supply. Install the new DC Power Supply by reversing the disassembly sequence.

After installation and power-up, it may be necessary to re-adjust the contrast on the front panel display (it may be too light or too dark)

9.8 Adjusting the 5 Volt DC Supply

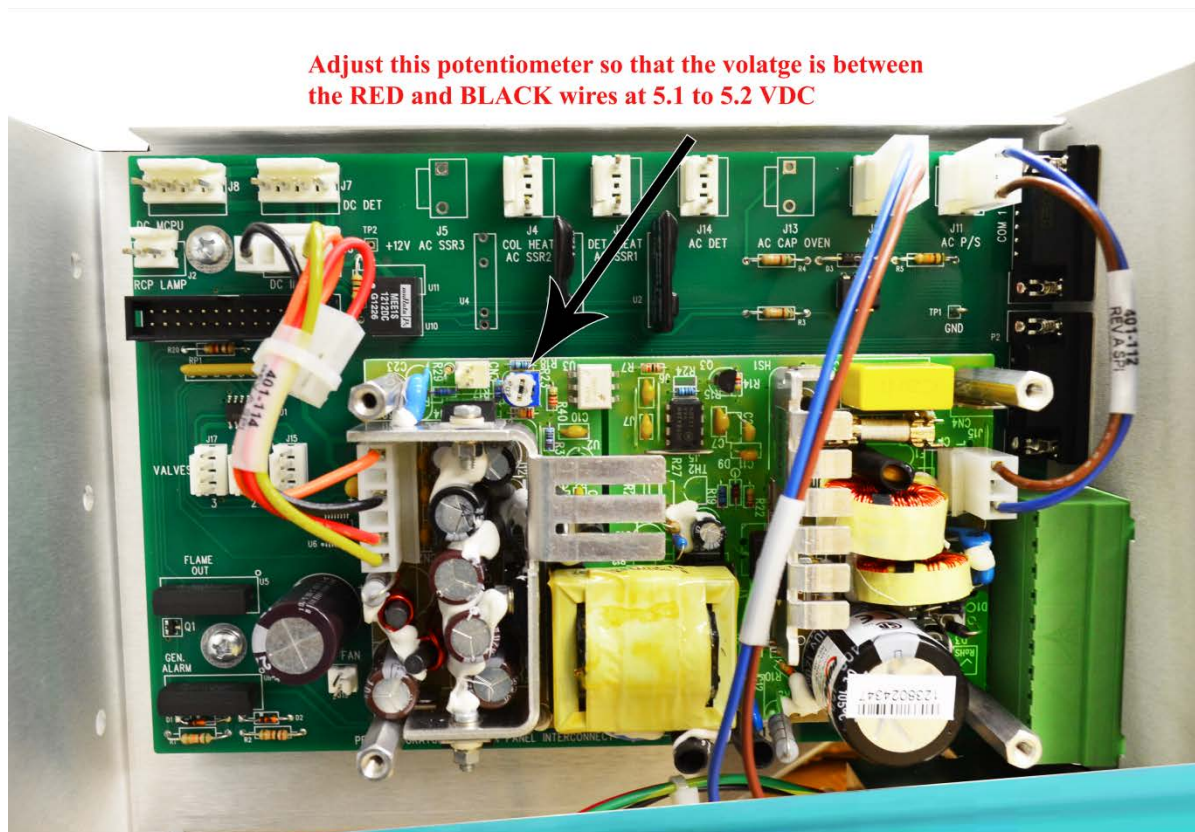

* Please Note: Adjustments is not applicable for newer analyzers (Rev L. and above). Please contact your local representative for details.

Figure 82. Adjusting the 5 Volt DC Supply

9.9 COM 1 Connections on Multiple Analyzers to a Single Computer

Equipment:

1 each - Computer: 1.2 GHz, 20 Gb Hard drive, 256 Mb RAM, CD-R drive, with at least 1 USB port:

(Dell Inspiron 8500 laptop or equivalent)

Operating System: Windows XP Home Edition, v2002 Service Pack 1

1 each - USB Hub: Inland "4 Port USB Hub", Self-powered, Stock #08302 (manufactured by Inland International – www.inland-products.com)

2 each - USB to DB9 (RS232) Adapter: SIIG, Inc "USB to Serial Adapter" Model # US2308, Part # JU-CB1S12 (manufactured by SIIG, Inc – www.siiq.com) *Comes with installation software*

2 each - DB9 Extension Cable: DB9M to DB9F, wired straight through, 10 feet long: GQ Cables "10' Mouse/Keyboard Extension, DB9M to DB9F" Part # M05-103

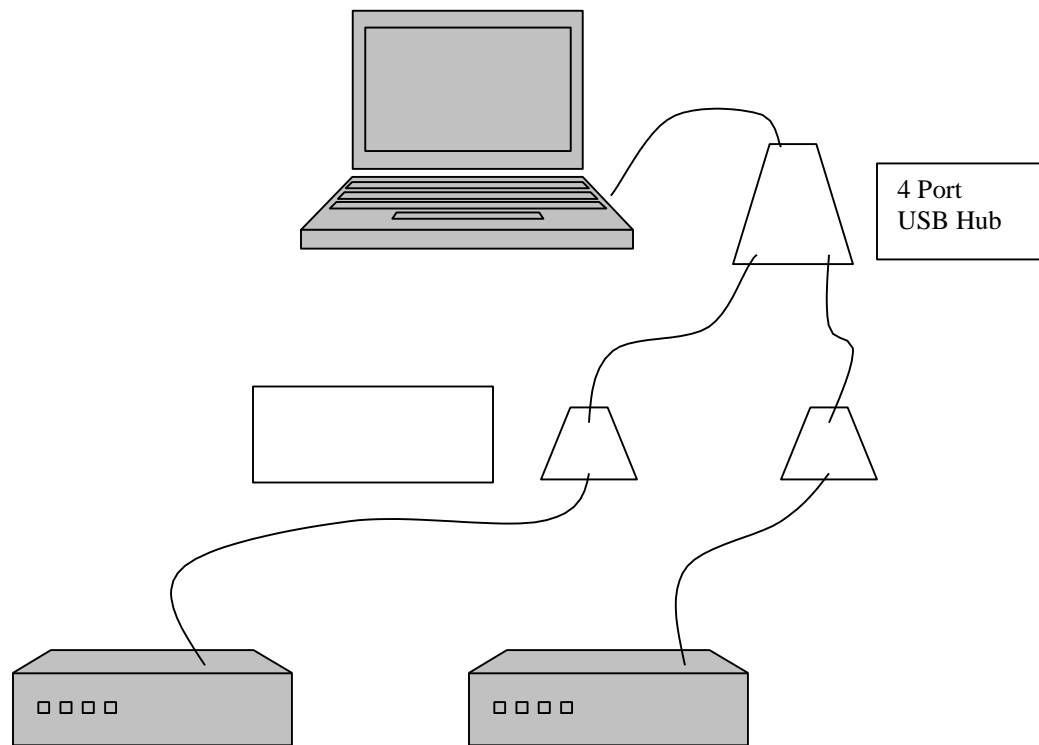


Figure 83. Multiple Analyzers Connected to a Single Computer For Using Viewer Software

Installation of Multiple Analyzers via USB:

- 1) Connect DB9M to DB9F extension cables to analyzers.
- 2) Power up analyzers – do not heat.
- 3) Set Port 1 Protocol on the analyzers to Viewer.

- 4) Insert 4 port USB hub into computer USB port – your computer should detect the connection and automatically install any required Windows drivers for USB hubs.
- 5) Connect the USB to DB9 adapters to the DB9 extension cables.
- 6) Insert first USB to DB9 adapter into 4 port USB hub – your computer should detect the connection and launch the hardware installation program.
- 7) Insert USB to DB9 adapter INSTALL CD in CD drive (it will launch automatically if auto play is enabled on your computer)
- 8) Follow installation instructions that come with USB to DB9 adapter (I confirm that it is okay to click the “Continue anyway” box)
- 9) Your computer may prompt to install the USB to DB9 adapter a second time – do not cancel, it is best to let the computer re-install the software again.
- 10) Create a new Desktop folder for each analyzer that you will connect with a USB to DB9 adapter. Be sure to give the folders a name that relates to the analyzer that you will connect to that adapter.
- 11) Copy the Viewer .exe files into each new folder. Open one of the new folders and select the Viewer.exe file
- 12) Double-click to launch a Viewer window.
- 13) Go to Windows' Device Manager (Settings>Control Panel>System>Hardware>Device Manager) to see what the port assignments on your computer are.
- 14) In the Device Manager, click on “Ports” to expand heading so you can see what COM port number your computer gave to the USB to DB9 adapter
- 15) Note which COM port number is now assigned to the new USB to Serial (RS232) Port (adapter) – it's COM4 on most computers.
- 16) Go back to the Viewer window
- 17) Under “File” on toolbar, select “Port Settings”
- 18) Enter COM number assigned to first USB to DB9 adapter (probably COM4), check that Baud Rate = 38400, Data Bits =8, Stop Bits = 1
- 19) Click OK box to save these settings.
- 20) Communications settings for that USB to DB9 adapter will be saved in the folder and will be used automatically every time you launch the Viewer software in that folder.
- 21) Leave the first Viewer window open.
- 22) Insert second USB to DB9 adapter into USB hub – your computer should detect the connection and launch the hardware installation program.
- 23) Insert USB to DB9 adapter INSTALL CD in CD drive (it will launch automatically if auto play is enabled on your computer)
- 24) Follow installation instructions that come with USB to DB9 adapter (I confirm that it is okay to click the “Continue anyway” box)
- 25) Your computer may prompt to install the USB to DB9 adapter a second time – do not cancel, it is best to let the computer re-install the software again.
- 26) Open the second new folder and select the Viewer.exe file in that folder.
- 27) Double-click to launch Viewer.exe (you might get a COM port error message window – simply click OK and proceed)
- 28) Go to Windows' Device Manager > Settings > Control Panel > System > Hardware > Device Manager) to see what the additional port assignment on your computer is.
- 29) In the Device Manager, click on “Ports” to expand heading so you can see what COM port number your computer gave to the second USB to DB9 adapter – it will be displayed directly under the note for the first USB to DB9 adapter.
- 30) Note which COM port number is now assigned to the second USB to Serial (RS232) Port (adapter) – it's COM5 on most computers.
- 31) Go back to the Viewer window

- 32) Under "File" on toolbar, select "Port Settings"
- 33) Enter COM number assigned to second USB to DB9 adapter (probably COM5), check that Baud Rate = 38400, Data Bits =8, Stop Bits = 1
- 34) Click OK box to save these settings.
- 35) Communications settings for that USB to DB9 adapter will be saved in the folder and will be used automatically every time you launch the Viewer software in that folder.
- 36) Viewer should automatically connect to the other analyzer.

At this point, you should have two PeakHost windows open, connected to the different analyzers.

If you run the analyzers, the chromatographic data for each analyzer will be automatically stored in the folder that matches where the Viewer.exe file is stored.