

*La Pierres*

The following is a user manual guideline for AngelTouch device by La Pierres

The AngelTouch is elegantly designed with ultrasonic multi-modulation-technology (MMT). Through its recommended use, it is your solution to more beautiful looking skin.

In addition to the advantages the AngelTouch features, you also have the option to choose a variety of La Pierres product lines to apply to your skin.

AngelTouch is also used and endorsed by skin care professionals worldwide.

## ABOUT ANGELTOUCH

AngelTouch is beautifully designed facial rejuvenating device equipped with La Pierres' patented Multi-Modulation-Technology (MMT). It enables deeper penetration and better absorption of the active ingredients in the skin care serum and cream to enhance their healing and rejuvenating powers. It delivers younger, healthier and better looking skin than using skin care products alone. It helps users realize a greater potential of the skin care products of their own choice.



Front



Side



Back

BASIC KIT

## KEY FEATURES (FRONT SIDE)

Medical grade smooth steel surface designed to glide on face with minimal effort

Direct dispensing avoid the need to apply serum or cream by finger tips and is more sanitary

Touch sensor control, with LED and sound indicators

Wireless inductive base charger




## KEY FEATURES (BACK SIDE)



## CONNECT FOR CHARGING

Connect wireless charger base to AC adaptor by plugging the smaller end of the USB cable into the inlet at the back of the wireless charger base with USB cable bigger end connecting to AC adaptor, and plug AC Adaptor into wall outlet

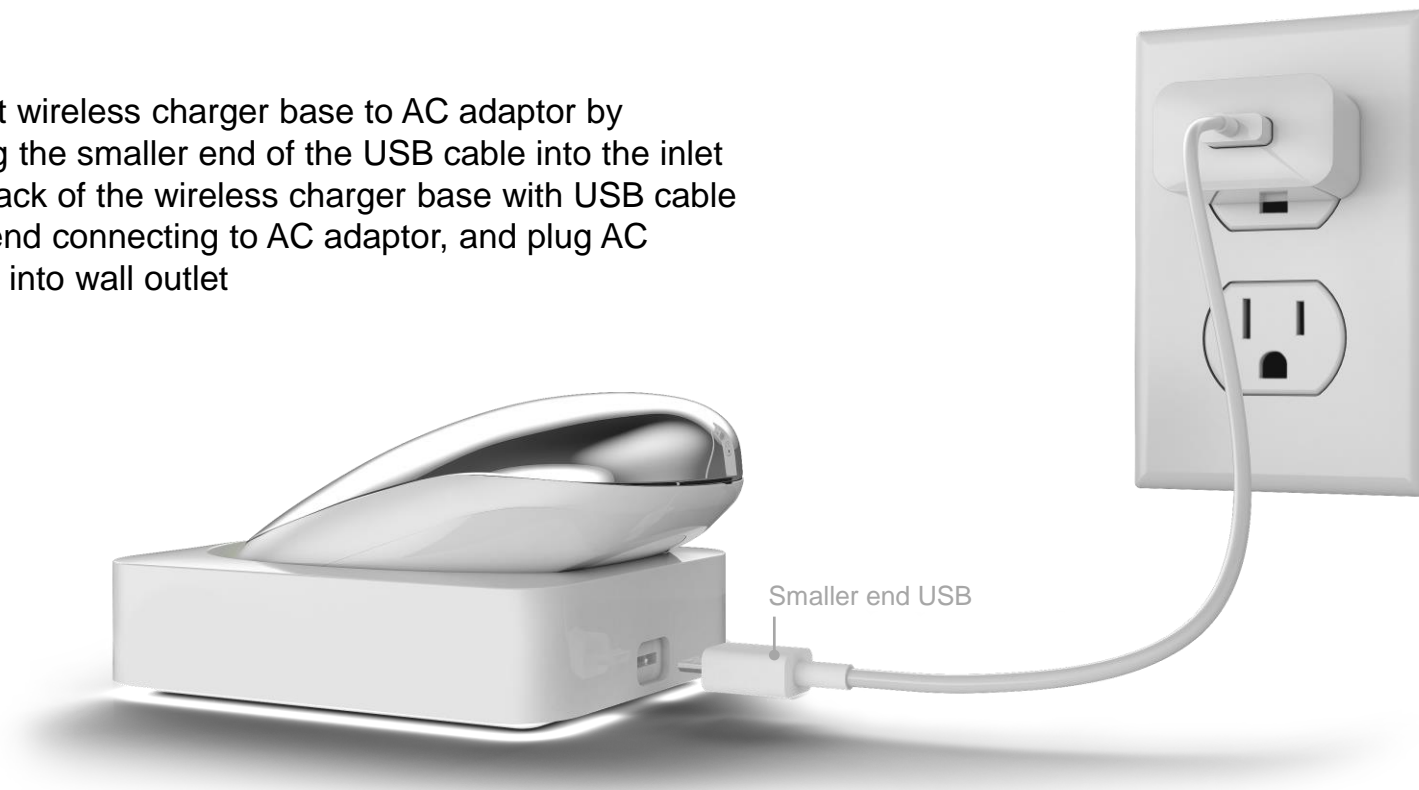

## AUTOMATIC AND EASY CHARGING


Place the AngelTouch device onto the wireless charger base as shown

Fully charged battery requires a first time charge of 8 hours.

*i.* When AngelTouch device is positioned on the wireless charger base, steady white light glows underneath.

*ii.* While charging, AngelTouch front logo will flash green light


When the green light on the AngelTouch device goes off, it is fully charged

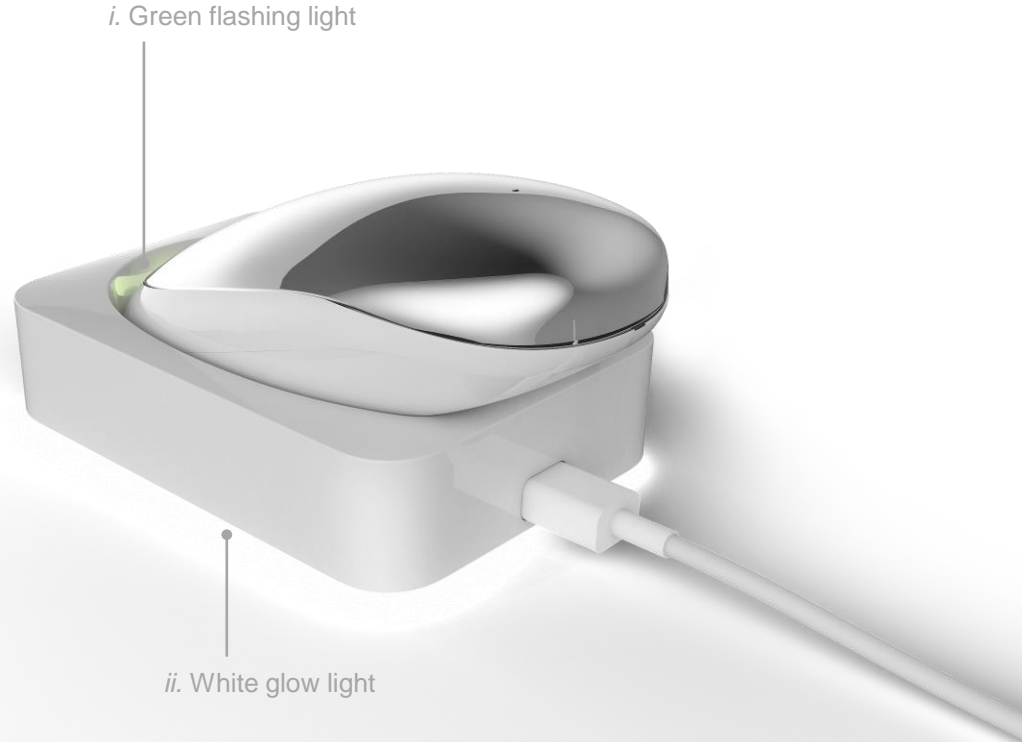



## 1 – 2 – 3 EASY INSTALL CARTRIDGE

- ① Push in the release button to remove the back cover
- ① Hold the cartridge by the side tabs and align the rail with the AngelTouch device groove. Holding the side ears position, push down until it stops
- ② Put back the cover at a forward angle and push the button end of the cover down until you hear “click”. It is now in lock position and ready for use

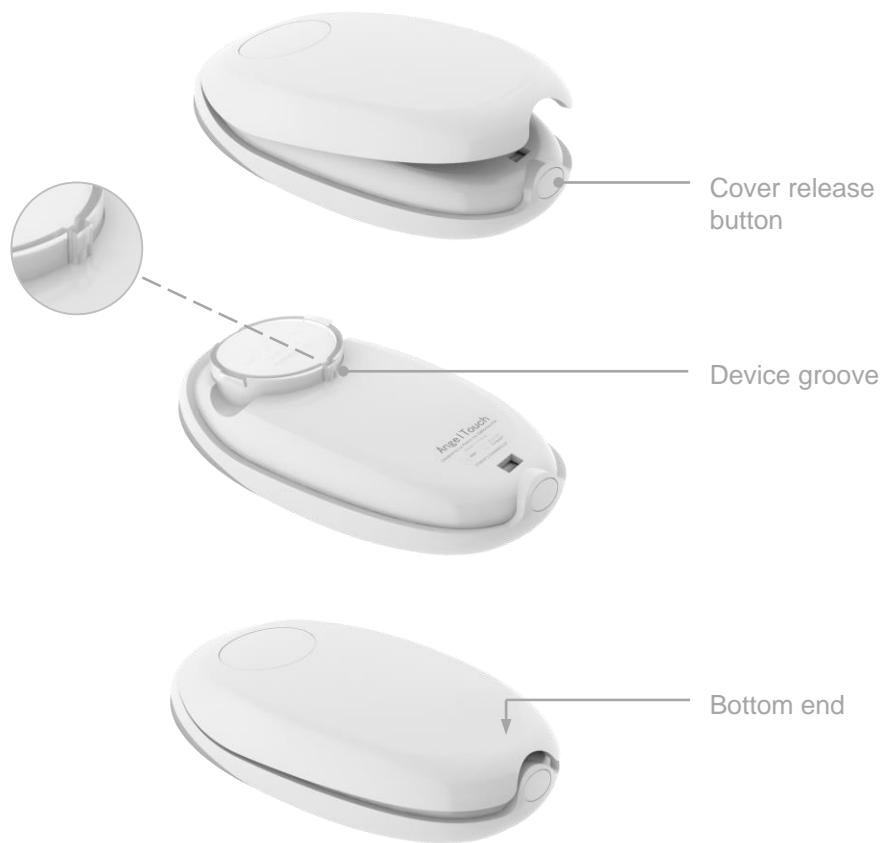

# POWERING ON & OFF

## TO TURN ON

Double tap on the logo located in the front of the AngelTouch device. The device will **beeps once** and the green light is on indicating that it is ready for use.

## TO TURN OFF

Double tap on the logo located in the front of the AngelTouch device. The device will **beeps twice** and green light goes off indicating it is powered off.

AngelTouch device is optimized to automatically turn off after 5 minutes of treatment for best results.

After 5 minutes, device will turn off, and beeps twice as an indication. Logo area green light goes off.



# ADVICE ON USING ANGELTOUCH

## FACIAL TREATMENT STEP

- ① Double tap logo to turn on
- ② Position metal surface of the device upon desired facial area
- ③ Make slow circular movement to spread serum/cream over desired facial area
- ④ Keep device steady over the desired facial area for 2~5 minutes
- ⑤ A slow rocking movement is effective in treating a larger area and enhancing absorption
- ⑥ After 5 minutes of usage, device will automatically shut down with 2 beeps as an indicator of finishing one treatment process

# RECOMMENDED APPLICATIONS

## DAILY REGIMEN AND SUGGESTED USE FOR ANGELFORMULA AND ANGELFORMULA CLEANSER

Direction for [AngelFormula Cleanser](#) before using AngelTouch device

*Dispense a moderate amount into palm of the hand and gently massage face (and neck) in circular motions*

*Remove with warm water and pat dry*

[AngelFormula 1](#) – All Day Anti-Aging Serum AM PM

For best results, apply twice per day, one in the morning and one at night

[AngelFormula 2](#) – Day Cream AM

Apply daily in the morning, layered with AngelFormula 1 for enhanced (anti-aging) treatment

[AngelFormula 3](#) – Night Serum PM

Apply daily at bedtime, layered

[Always apply AngelFormula creams or serums on clean skin and enjoy the results!](#)

## DEVICE GUIDELINES

- Metal plate warming up is normal during treatment. This is due to acoustic wave being activated
- After usage the skin may continue to enjoy the invigorating sensation from the use of the device and applied skin care product. Skin may also experience renewed hydration and smoother texture.
- A fully charged device can provide up to 90 minutes treatment time  
**Recommended:** Maximum treatment time on any specific spot of facial area in a single treatment is **5 minutes**
- One facial treatment should not exceed **15 minutes**;  
2 treatments per day is preferred and maximum usage per day is **30 minutes**

**More usage than what is recommended may cause discomfort on sensitive skin**

## DEVICE GUIDELINES CONTINUED

- After treatment, the device metal surface may be cleaned with a damp cloth.  
DO NOT WASH
- Before charging the device, the cartridge should always be removed from device. Put cartridge back into the storage box provided
- Important: After any one day of use, the dispensing outlet of the device should be thoroughly cleaned with the cotton swabs provided. This will remove any residue inside and outside of the device, avoiding a potential clog and maintains optimal performance. The dispenser should be cleaned with alcohol damped swabs at least once per week
- Device should always be re-charged to full power after extended idle time before the next usage

## SKIN CARE GUIDELINES

- Always use AngelFormula or La Pierres recommended products with AngelTouch
- La Pierres skin care contains pure, potent botanicals combined with scientifically advanced ingredients.
- Application of our effective formulations with the AngelTouch Device promotes enhanced, healthy radiant skin
- AngelFormula have been designed specifically for the AngelTouch device. The AngelFormula will achieve the best results with our pure natural high quality serums
- Majority of non-La Pierres skin care products contain both active and inactive ingredients
- AngelTouch is effective by enhance absorption of the active ingredients. Inactive additives may form paste-like or oil-like layer, which is not a defect of the device

[Check out our other line of products for better results with AngelTouch device](#)

# LIGHT AND SOUND INDICATIONS

## AFTER DOUBLE-TAPPING LOGO AT THE BOTTOM OF DEVICE...

- 1 beep + steady green light: *Device is turned on correctly*
- 2 beeps + green light turn off: *Device is turned off correctly*
- 3 red flashes: *Low battery: Charge device now*
- 10 red flashes: *Wait for device to cool down*

## WHEN YOU ARE USING DEVICE...

- 2 beeps + green light turn off: *Auto turn-off after 5 minutes usage*
- Repeated 1 beep: *Battery near empty, prepare to end treatment*
- 3 red flashes + 2 beeps: *Auto turn-off due to low battery. Charge device now*
- 10 red flashes + 2 beeps: *Auto turn-off. Wait for device to cool down*

## WHEN DEVICE IS CHARGING...

- Slow glowing green light: *Normal charging*
- Glowing green light turns off: *Battery is fully charged*



# SPECIFICATIONS

AC adaptor: 100~240V, 50~60Hz

Charge time: 4 hours (from empty to full)

Use time: 18 times of 5-minutes treatment <sup>(a)</sup>

Non-use standby time: 1 week from full charge <sup>(a)</sup>

Auto turnoff: 5 minutes

Operation Temperature: -10~40°C

(a) Maximum time for new device.

Actual time may vary due to environment difference and battery aging over lifetime

## REGISTER YOUR PRODUCT

Register your AngelTouch device so that we can better serve you. You can sign up for La Pierres email updates and gain access to private invitations for promotions and new product information.

# CAUTION

## WARNINGS: FOR ANGELTOUCH AND INCLUDED ACCESSORIES

DO NOT submerge in water

DO NOT allow water to get into the device or the charger base

DO NOT expose to fire

DO NOT drop on hard surface

DO NOT expose to temperature higher than 104°F (40°C) or lower than 14°F (-10°C)

KEEP OUT of the reach of children

DO NOT use in ways other than described in this manual

MUST follow cleaning procedures for the device as described in this manual

BREACHING ANY OF THE ABOVE TERMS WILL VOID THE WARRANTY OF THE ANGELTOUCH AND INCLUDED ACCESSORIES.

# WARRANTY

## LA PIERRES ONE (1) YEAR LIMITED WARRANTY

LA PIERRES WARRANTY OBLIGATIONS FOR THIS ELECTRONIC PRODUCT ARE LIMITED TO THE TERMS SET FORTH BELOW: LA PIERRES INC. ("LA PIERRES") WARRANTS THIS LA PIERRES-BRANDED ELECTRONIC PRODUCT AGAINST DEFECTS IN MATERIALS AND WORKMANSHIP UNDER NORMAL USE FOR A PERIOD OF ONE (1) YEAR FROM THE DATE OF RETAIL PURCHASE BY THE ORIGINAL END-USER PURCHASER ("WARRANTY PERIOD"). IF A ELECTRONIC DEFECT ARISES AND A VALID CLAIM IS RECEIVED WITHIN THE WARRANTY PERIOD, AT ITS OPTION AND TO THE EXTENT PERMITTED BY LAW, LA PIERRES WILL EITHER (1) REPAIR THE PRODUCT DEFECT AT NO CHARGE, USING NEW OR REFURBISHED REPLACEMENT PARTS, (2) EXCHANGE THE PRODUCT WITH A PRODUCT THAT IS NEW OR WHICH HAS BEEN MANUFACTURED FROM NEW OR SERVICEABLE USED PARTS AND IS AT LEAST FUNCTIONALLY EQUIVALENT TO THE ORIGINAL PRODUCT. LA PIERRES MAY REQUEST THAT YOU REPLACE DEFECTIVE PARTS WITH NEW OR REFURBISHED USER-INSTALLABLE PARTS THAT LA PIERRES PROVIDES IN FULFILLMENT OF ITS WARRANTY OBLIGATION. A REPLACEMENT PRODUCT OR PART, INCLUDING A USER-INSTALLABLE PART THAT HAS BEEN INSTALLED IN ACCORDANCE WITH INSTRUCTIONS PROVIDED BY LA PIERRES, ASSUMES THE REMAINING WARRANTY OF THE ORIGINAL PRODUCT OR NINETY (90) DAYS FROM THE DATE OF REPLACEMENT OR REPAIR, WHICHEVER PROVIDES LONGER COVERAGE FOR YOU. WHEN A PRODUCT OR PART IS EXCHANGED, ANY REPLACEMENT ITEM BECOMES YOUR PROPERTY AND THE REPLACED ITEM BECOMES LA PIERRES PROPERTY. PARTS PROVIDED BY LA PIERRES IN FULFILLMENT OF ITS WARRANTY OBLIGATION MUST BE USED IN PRODUCTS FOR WHICH WARRANTY SERVICE IS CLAIMED. WHEN A REFUND IS GIVEN, THE PRODUCT FOR WHICH THE REFUND IS PROVIDED MUST BE RETURNED TO LA PIERRES AND BECOMES LA PIERRES PROPERTY.

# WARRANTY CONTINUED

## EXCLUSIONS AND LIMITATIONS

THIS LIMITED WARRANTY APPLIES ONLY TO ELECTRONIC PRODUCTS MANUFACTURED BY OR FOR LA PIERRES THAT CAN BE IDENTIFIED BY THE “LA PIERRES” OR “ANGELTOUCH” TRADEMARK, TRADE NAME, OR LOGO AFFIXED TO THEM. THE LIMITED WARRANTY DOES NOT APPLY TO ANY NON-LA PIERRES ELECTRONIC PRODUCTS OR ANY NON-ELECTRONIC PRODUCTS, EVEN IF PACKAGED OR SOLD BY LA PIERRES. MANUFACTURERS, SUPPLIERS, OR PUBLISHERS, OTHER THAN LA PIERRES, MAY PROVIDE THEIR OWN WARRANTIES TO THE END USER

PURCHASER, BUT LA PIERRES, IN SO FAR AS PERMITTED BY LAW, PROVIDES THEIR PRODUCTS “AS IS”. NON-ELECTRONIC PRODUCTS DISTRIBUTED BY LA PIERRES WITH OR WITHOUT THE LA PIERRES BRAND NAME IS NOT COVERED UNDER THIS LIMITED WARRANTY. LA PIERRES DOES NOT WARRANT THAT THE OPERATION OF THE PRODUCT WILL BE UNINTERRUPTED OR ERROR-FREE. LA PIERRES IS NOT RESPONSIBLE FOR DAMAGE ARISING FROM FAILURE TO FOLLOW INSTRUCTIONS RELATING TO THE PRODUCT’S USE. THIS WARRANTY DOES NOT APPLY: (A) TO DAMAGE CAUSED BY USE WITH NON-LA PIERRES PRODUCTS; (B) TO DAMAGE CAUSED BY ACCIDENT, ABUSE, MISUSE, FLOOD, FIRE, EARTHQUAKE OR OTHER EXTERNAL CAUSES; (C) TO DAMAGE CAUSED BY OPERATING THE PRODUCT OUTSIDE THE PERMITTED OR INTENDED USES DESCRIBED BY LA PIERRES; (D) TO DAMAGE CAUSED BY SERVICE (INCLUDING UPGRADES AND EXPANSIONS) PERFORMED BY ANYONE WHO IS NOT A REPRESENTATIVE OF LA PIERRES OR A LA PIERRES AUTHORIZED SERVICE PROVIDER;

(E) TO A PRODUCT OR PART THAT HAS BEEN MODIFIED TO ALTER FUNCTIONALITY OR CAPABILITY WITHOUT THE WRITTEN PERMISSION OF LA PIERRES; (F ) TO CONSUMABLE PARTS, SUCH AS DISPENSERS, UNLESS DAMAGE HAS OCCURRED DUE TO A DEFECT IN MATERIALS OR WORKMANSHIP; (G) TO COSMETIC DAMAGE, INCLUDING BUT NOT LIMITED TO SCRATCHES, DENTS AND

## WARRANTY CONTINUED

OR (H) IF ANY LA PIERRES SERIAL NUMBER HAS BEEN REMOVED OR DEFACED. THIS WARRANTY GIVES YOU SPECIFIC LEGAL RIGHTS, AND YOU MAY ALSO HAVE OTHER RIGHTS WHICH VARY FROM STATE TO STATE. LA PIERRES' RESPONSIBILITY FOR ELECTRONIC DEFECTS IS LIMITED TO REPAIR OR REPLACEMENT SERVICE AS DETERMINED BY LA PIERRES IN ITS SOLE DISCRETION. ALL EXPRESS AND IMPLIED WARRANTIES, INCLUDING BUT NOT LIMITED TO ANY IMPLIED WARRANTIES AND CONDITIONS OF MERCHANTABILITY AND FITNESS FOR A

PARTICULAR PURPOSE, ARE LIMITED IN TIME TO THE TERM OF THIS LIMITED WARRANTY. NO WARRANTIES, WHETHER EXPRESS OR IMPLIED, WILL APPLY AFTER THE EXPIRATION OF THE LIMITED WARRANTY PERIOD. SOME STATES DO NOT ALLOW LIMITATIONS ON HOW LONG AN IMPLIED WARRANTY LASTS, SO THIS LIMITATION MAY NOT APPLY TO YOU. NO LA PIERRES RESELLER, AGENT, OR EMPLOYEE IS AUTHORIZED TO MAKE ANY MODIFICATION, EXTENSION, OR ADDITION TO THIS WARRANTY. IF ANY TERM IS HELD TO BE ILLEGAL OR UNENFORCEABLE, THE LEGALITY OR ENFORCEABILITY OF THE REMAINING TERMS SHALL NOT BE AFFECTED OR IMPAIRED. EXCEPT AS PROVIDED IN THIS WARRANTY AND TO THE EXTENT PERMITTED BY LAW, LA PIERRES IS NOT RESPONSIBLE FOR DIRECT, SPECIAL, INCIDENTAL OR CONSEQUENTIAL DAMAGES RESULTING FROM ANY BREACH OF WARRANTY OR CONDITION, OR UNDER ANY OTHER LEGAL THEORY, INCLUDING BUT NOT LIMITED TO LOSS OF USE; LOSS OF REVENUE; LOSS OF ACTUAL OR ANTICIPATED PROFITS (INCLUDING LOSS OF PROFITS ON CONTRACTS); LOSS OF THE USE OF MONEY; LOSS OF ANTICIPATED SAVINGS; LOSS OF BUSINESS; LOSS OF

OPPORTUNITY; LOSS OF GOODWILL; LOSS OF REPUTATION; LOSS OF, DAMAGE TO OR CORRUPTION OF DATA; OR ANY INDIRECT OR CONSEQUENTIAL LOSS OR DAMAGE HOWSOEVER CAUSED INCLUDING THE REPLACEMENT OF EQUIPMENT AND PROPERTY, ANY COSTS OF RECOVERING, PROGRAMMING, OR REPRODUCING ANY USER DATA STORED OR USED WITH LA PIERRES PRODUCTS AND ANY FAILURE TO MAINTAIN THE CONFIDENTIALITY OF DATA STORED ON THE PRODUCT. THE FOREGOING LIMITATION SHALL NOT

## WARRANTY CONTINUED

PERSONAL INJURY CLAIMS, OR ANY STATUTORY LIABILITY FOR INTENTIONAL AND GROSS NEGLIGENT ACTS AND/OR OMISSIONS. SOME STATES DO NOT ALLOW THE EXCLUSION OR LIMITATION OF INCIDENTAL OR CONSEQUENTIAL DAMAGES, SO THE ABOVE LIMITATION MAY NOT APPLY TO YOU.

### OBTAINING WARRANTY SERVICE

PLEASE ACCESS AND REVIEW THE ONLINE HELP RESOURCES REFERRED TO IN THE DOCUMENTATION ACCOMPANYING THIS ELECTRONIC PRODUCT BEFORE REQUESTING WARRANTY SERVICE. IF THE PRODUCT IS STILL NOT FUNCTIONING PROPERLY AFTER MAKING USE OF THESE RESOURCES, PLEASE CONTACT THE LA PIERRES REPRESENTATIVES OR, IF APPLICABLE, AN LA PIERRES LOCAL OFFICE OR LA PIERRES AUTHORIZED SERVICE PROVIDER LOCATED USING THE INFORMATION PROVIDED IN THE DOCUMENTATION. WHEN CONTACTING LA PIERRES VIA TELEPHONE, OTHER CHARGES MAY APPLY BECOMES THE PROPERTY OF LA PIERRES AND YOU AGREE TO FOLLOW INSTRUCTIONS, INCLUDING, IF REQUIRED, ARRANGING THE RETURN OF ORIGINAL PRODUCT OR PART TO LA PIERRES IN A TIMELY MANNER. WHEN PROVIDING DIY SERVICE REQUIRING THE RETURN OF THE ORIGINAL PRODUCT OR PART, LA PIERRES MAY REQUIRE A CREDIT CARD AUTHORIZATION AS SECURITY FOR THE RETAIL PRICE SERVICE OPTIONS, PARTS AVAILABILITY AND RESPONSE TIMES MAY VARY. SERVICE OPTIONS ARE SUBJECT TO CHANGE AT ANY TIME. IN ACCORDANCE WITH APPLICABLE LAW, LA PIERRES MAY REQUIRE THAT YOU FURNISH PROOF OF PURCHASE DETAILS AND/OR COMPLY WITH REGISTRATION REQUIREMENTS BEFORE RECEIVING WARRANTY SERVICE. PLEASE REFER TO THE ACCOMPANYING DOCUMENTATION FOR MORE DETAILS ON THIS AND OTHER MATTERS ON OBTAINING WARRANTY

## WARRANTY CONTINUED

SERVICE. LA PIERRES WILL MAINTAIN AND USE CUSTOMER INFORMATION IN ACCORDANCE WITH THE LA PIERRES CUSTOMER PRIVACY POLICY AT [WWW.LAPIERRES.COM/SUPPORT/PRIVACY](http://WWW.LAPIERRES.COM/SUPPORT/PRIVACY).

# DISCLAIMER

LA PIERRES INC. CONTINUOUSLY STRIVES TO ENHANCE ITS PRODUCTS EFFECTIVENESS IN IMPROVING USERS' FACIAL SKIN CONDITION AND TO MEET EXPECTATIONS OF ALL OUR CUSTOMERS. HOWEVER, AS HUMAN SKIN CONDITION AND SKIN SENSITIVITY ARE UNPREDICTABLE AND SHOW LARGE VARIATIONS AMONG DIFFERENT HUMAN GROUPS, IT IS EXTREMELY DIFFICULT TO ACHIEVE THE SAME KIND OF SKIN IMPROVEMENT RESULT FOR ALL OUR CUSTOMERS. IT IS IMPOSSIBLE FOR LA PIERRES TO ASSESS THE FULL RANGE OF SKIN CONDITIONS AND SKIN SENSITIVITY LEVELS OF OUR USERS. THUS, LA PIERRES IS NOT LIABLE FOR ANY CLAIM THAT THE RESULT OF SKIN TREATMENT IS NOT PRODUCING CLAIMED EFFECTS. LA PIERRES IS NOT LIABLE FOR THE RARE CASES OF SKIN CONDITION DEGRADE OR DAMAGE AFTER USING LA PIERRES PRODUCTS, OR ANY BODILY OR PROPERTY DAMAGES OR MONETARY LOSSES INCURRED. IT IS THE RESPONSIBILITY OF THE USERS TO TEST OUT THE LA PIERRES PRODUCTS COMPATIBILITY AND SAFETY WITH THEIR OWN SKIN CONDITIONS BEFORE USAGE. ADDITIONALLY, DUE TO THE AMOUNT OF VARIATION IN THE SKIN CARE PRODUCTS (SERUM, CREAM, ETC.), LA PIERRES DOES NOT GUARANTEE THAT ANY SKIN CARE PRODUCT THAT USERS CHOOSE WOULD BE COMPATIBLE WITH LA PIERRES PRODUCTS AND PRODUCE RESULTS THAT USERS DESIRE. LA PIERRES IS NOT LIABLE FOR ANY LOSS INCURRED BY NON-COMPATIBLE SKIN CARE PRODUCTS. LA PIERRES IS ALSO NOT LIABLE FOR BODILY OR PROPERTY DAMAGES CAUSED BY ANY IMPROPER USAGE OF LA PIERRES PRODUCTS, OR BY USING LA PIERRES PRODUCTS WITH ANY NON-LA PIERRES PRODUCTS TOGETHER OR SEPARATELY.



For A Touch of Beauty

designed by

*La Pierres*

©2015 La Pierres, Inc. All rights reserved.

La Pierres, AngelTouch, AngelFormula and La Pierres logo are trademarks of

La Pierres Inc., registered in the U.S. and other countries.

Designed by La Pierres in San Jose, California USA. Ver.2015-09