




68 Route 125
Kingston, New Hampshire 03848
U.S.A.
T. (603) 347-3005 F. (603) 642-9291
www.infinitetherapeutics.com

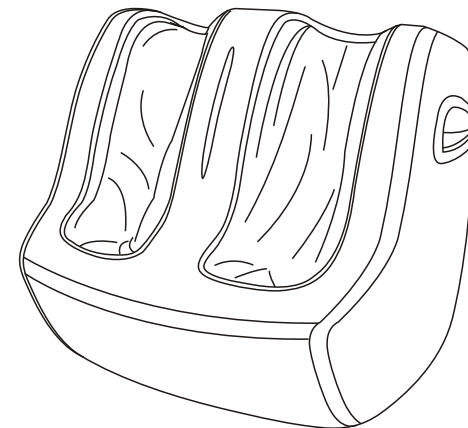

Leg & Foot Massager

RT1880
USER MANUAL

Thank you for your purchase of this product. Before use, please read this USER MANUAL carefully, and pay special attention to the safety precautions to ensure proper operation and use.

Please keep this USER MANUAL properly for further reference.

CONTENTS

SAFETY PRECAUTIONS-----	1
NAME AND FUNCTION OF EACH PART-----	3
External Structure-----	3
Internal Structure-----	4
Functions on Control Panel-----	5
USE METHOD-----	6
Preparation before Use-----	6
Operation Method-----	6
Use method-----	7
After Use-----	7
CLEANING AND MAINTENANCE-----	8
Maintenance-----	8
Kind reminders-----	8
TROUBLE SHOOTING-----	9
TECHNICAL SPECIFICATIONS-----	9

TROUBLESHOOTING

No.	Symptom	Possible Cause	Cure
1	There is a clattering noise in the movement.	This is perfectly normal caused by the working of product.	It is normal and no need for treatment.
2	The Foot Massager can not be started up.	The power plug is loosened.	Insert the plug into socket firmly.
		The power switch is not turned on.	Turn on the power switch.
		Fuse is burnt	Change fuse
3	Stop during use	Default working time is up	Reuse in 30 minutes later

If there are some unusual symptoms after having checked the above items, please contact your local dealer to ask for an authorized service engineer.

TECHNICAL SPECIFICATIONS

Product Name:	INFINITY DELUXE LEG & FOOT MASSAGER
Model No.:	IT-1880
Rated Voltage:	220-240 V~
Rated Frequency:	50/60Hz
Rated Power:	145 W
Safety Design:	Class I

► CLEANING AND MAINTENANCE

Maintenance

1. Cleaning

- Before cleaning, make sure the power supply is turned off and the plug has been pulled out from the socket.
- Use a soft cloth to wipe the surface of the product and controller.
- Use a piece of cloth that has been dipped with cleaning agents to wipe the main unit of massage chair, and then use a piece of dry and soft cloth to wipe it.

⚠ Caution

- Do not pull the power plug by damp hands.
- Do not use water washing product, do not spray the water or other liquids to the interior of the product.
- Do not use corrosive cleaning agents like benzene or thinner to clean the massage chair, otherwise, the color of upholstery will fade or the massage chair will crack or be eroded.

2. Maintenance

- Make sure the power supply is turned off and the plug has been pulled out from the socket.
- Sweep off the dust. Keep the massage chair clean and away from damp places.
- If the massage chair is not to be used for a long period of time, please cover the massage chair with a dustproof cover (self provided).

⚠ Caution

- Do not store the product under direct sunlight or in places with high temperature. This may cause the color of the upholstery to fade.

Kindly reminding

Please check if there is any below symptoms when the product is in daily use:

- The burned and scorched flavor
- Touch the wire, sometimes power is on, sometimes power is off
- The power cord heats.
- Other abnormalities.





⚠ Caution

- If any of the above mentioned abnormalities exist, please stop use to avoid accident happens.
- Turn off the power supply and pull the plug out, contact the local dealer or contact us.
- Do not disassemble the product and repair.




► SAFETY PRECAUTIONS

Safety Precautions

- Before using, please read this “SAFETY PRECAUTIONS” carefully to ensure proper operation.
- The following “SAFETY PRECAUTIONS” is of great help to your correctly use of this product.
- To clearly indicate the degrees of danger and damage, the classes of accidents caused by incorrect use are classified into “WARNING” and “CAUTION”, instructions marked as followings should be strictly followed:



	WARNING	Improper operation may cause injury to users.
	CAUTION	Improper operation may cause users injury or product damage.
	Forbidden. Specific instruction will be indicated in the diagram drawings.(e.g, the left drawings means disassembly is not allowed.)	
	Compulsory Specific instruction will be indicated in the diagram drawings.(e.g, the left drawings means to disconnect the power plug.)	

⚠ WARNING

 COMPULSARY	<p>We recommend that those with any of the following conditions consult a doctor before using.</p> <div><div><p>1. Those who are using electronic medical instrument imbedded in the body, such as cardiac pacemaker etc.</p><p>2. Those who are being treated by doctors, especially those who feel uncomfortable.</p><p>3. Those who are malignant tumor patients, heart disease patients, acute patients.</p></div><div><p>4. Women who are pregnant or menstruating.</p><p>5. Those with osteoporosis or fracture spine.</p><p>6. Those with skin disease or those whose skins are injured.</p><p>7. Those whose body temperature is over 38 degrees (Febrile Phase).</p></div></div>
 PROHIBITED	<p>Supervision is compulsory when the product is used by intellectual disabilities children.</p>
	<p>Do not allow children to play with product.</p>
	<p>Do not allow more than one person to use this product at a time.</p>
	<p>Do not damage, process, excessively bend, pull hard, twist or knot power cord.</p>
	<p>Do not use the product when the power cord is damaged.</p>
 NO DISASSEMBLY	<p>Do not allow the minors to use the product.</p>
	<p>Anyone else, except maintenance personnel specified by our company, should not disassemble, maintain, or remodel this product. If any problem, please contact local dealer immediately.</p>

SAFETY PRECAUTIONS

CAUTION

 COMPULSARY	The default working time of the product is 10 minutes
	Do not put heavy objects on the product
	Please do not cover the product with blanket or pillow when it is used, thus to avoid user from being burn or causing fire disaster.
	Please do not stand on the product to avoid any failure
	Do not operate the controller with wet hands.
	Do not put the product outdoors.
	It is forbidden to dip the whole product in water or other liquid when cleaning.
	If the product is moved from a low temperature place to a warmer place, it is advised to use it one hour later to avoid function failure caused from vapor condensing.
	Do not use this product together with other therapeutic equipments.
	Do not use this product for medical treatment.
	People who are insensitive to heat must be accompanied by normal person when using this product.
	Please stop using this product immediately when you find any abnormal situations (such as smoke or peculiar smell) during use.
	Please use the power supply with rated voltage only
	Please check the product before using, if damage found, please contact the dealer for repairing immediately
	Do not use this product in damp environment to avoid electric shock.
 PROHIBITED	During use, if you feel abnormal, please stop using this product immediately, and then consult your doctor.
	After use, please cut off power supply.
	Do not use the product within half an hour after meal
	Do not sleep during the product working.
	Do not use the product after drunk or any abnormality of your body.
	Do not wear watch, bangle and not put cell phone, keys and other valuables in your jacket when use the product.
	Do not use petrol, caustic agents and spray pesticide for cleaning.
	Do not put the product near combustible.
	Do not smoke when using the product to avoid fire occurrence

USE METHOD

Use method

Step 1:

- ▶ Connect the plug to rated power supply (220V~).



Caution

- ▶ Before connecting the plug to power supply, please check and make sure the power cord undamaged.
- ▶ Please do not insert or remove the plug with wet hands.

Step 2

- ▶ Turn on the power switch on back of the product, the LED indicator on the control panel will keep flashing.

Step 3

- ▶ Choose a chair with proper height and seat yourself, then put your feet into the product.



Caution

- ▶ Before use, please take valuable decorations off from feet to avoid damage or injury.
- ▶ Before use, please keep kids away from the product.

Step 4

- ▶ Choose the desirable massage program by using the control panel

After use

- ▶ Turn off the power switch on the back of this product
- ▶ Remove the plug from the power supply socket


► USE METHOD

Preparation before use

- Please put the product on flat floor.
- Keep this product 1 meter away from audio-visual equipment such as TV set, radio, etc to avoid signal interference.


⚠ Caution:

- When electrical failure or sudden stop working happens, please keep calm and then to open the rotary switch on the back of the product and rotary the switch clockwise or anticlockwise to release the splints. See images blow.
- After release, please close the cover onto the rotary switch, otherwise the product can not work.
- During using, the product will stop working if the switch cover is opened; it can only continue to use after the cover is closed.
- The rotary switch is working at the same time as product.

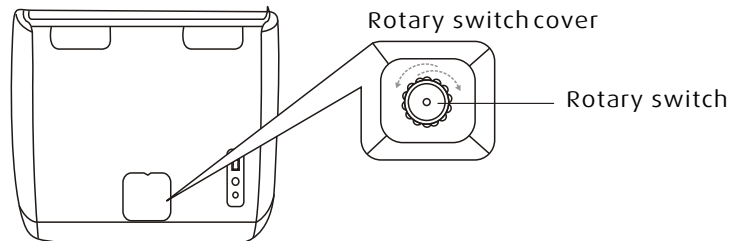

Moving method

- Grasp the handle on both side of the product to move the product. See image


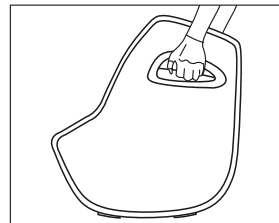


Caution

- Before moving the product, please switch off the power and remove the plug from the power socket.
- The floor may be scratched by the product, so it is advised to put a carpet under the bottom to protect the floor and reduce noise.


► NAME AND FUNCTION OF EACH PART

External structure

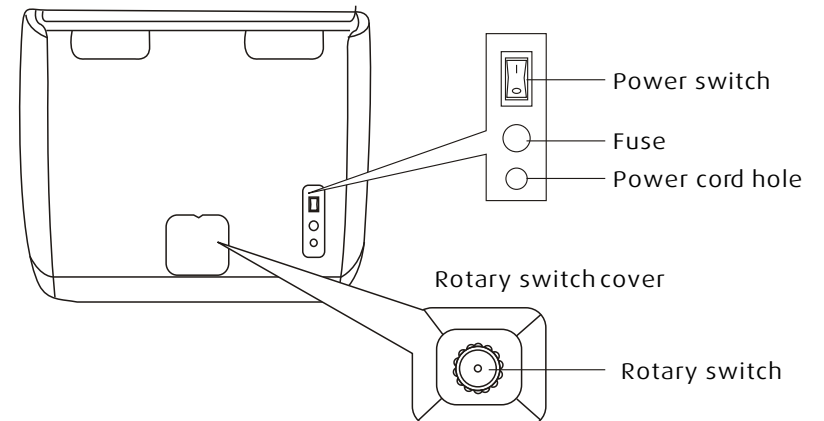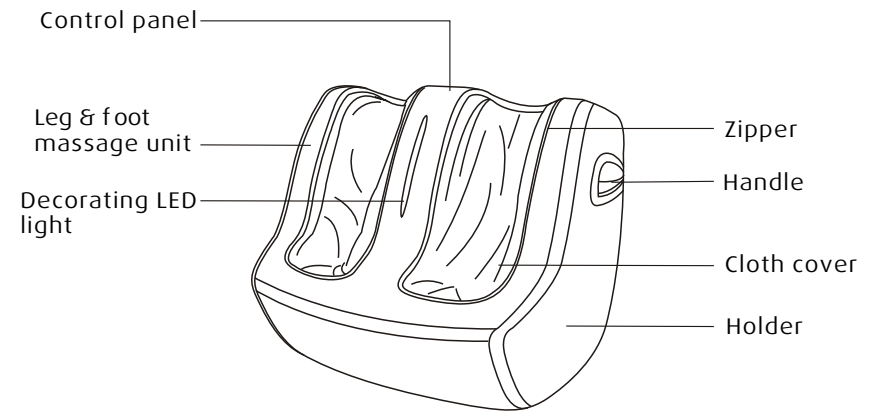


► NAME AND FUNCTION OF EACH PART

Internal structure

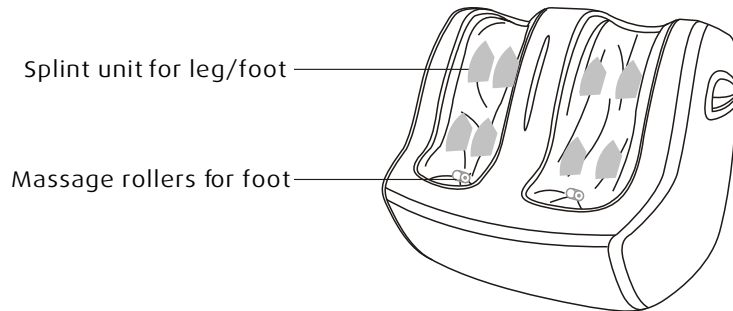

Control panel

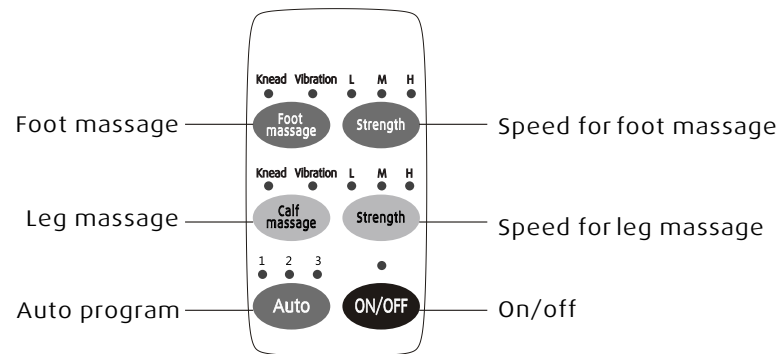

⚠ Caution

- It is not allowed to use the product by standing on it or put heavy objects on the product.
- In any case, please do not sit on the product.
- In case of electrical failure during using, please shut off the power switch and remove the plug from power supply.
- Please do not put hands or fingers into the gaps to avoid any injury
- Please keep kids and pets away from the product.
- Please do not use the product by no-load to prevent the cloth covers from being worn.

► NAME AND FUNCTION OF EACH PART

Functions description on control panel

- Foot massage key:**
Under working mode, press this key to enter to kneading/vibration massage function, press twice to switch to kneading function, press the key third time to start vibration function, press the key fourth time to stop kneading/vibration function, you can choose varies massage mode by successively cycle, and the corresponding LED indicator will be turned on or off.
- Speed key for foot massage:**
Under working mode, press foot massage key, the default speed is medium level; press this key to choose high, medium, low level successively, the corresponding LED indicator will be turned on or off.
- Leg massage key:**
Under working mode, press this key to enter to kneading/vibration massage function, press twice to switch to kneading function, press the key third time to start vibration function, press the key fourth time to stop kneading/vibration function, you can choose varies massage mode by successively cycle, and the corresponding LED indicator will be turned on or off.
- Speed key for leg massage:**
Under working mode, press leg massage key, the default speed is medium level; press this key to choose high, medium, low level successively, the corresponding LED indicator will be turned on or off.
- Auto key:**
Under standby mode, press this key to enter to auto massage function, auto 1, 2, 3 massage functions are available for choices, the corresponding LED indicator will be turned on or off.
- Power switch key:**
When powered on, the switch indicator will keep flashing, press this key for 1 second, the indicator will turn to light mode. The default massage program is auto 1 when the product starts working; press this key twice to stop working. (The default working time is 10 minutes)

⚠ Caution

- Please do not splash water or other liquids onto the control panel.
- Please do not keep pressing the keys with too much force to avoid any failure of program.
- The LED light is only flashing under product working mode.