UC-324THW/UC-324TMW/ UC-324THX Precision Scale Users Manual

English | French | Spanish

Table of Content

Before Use

- Precautions E-1
- Welcome E-2
- Scale Components..... E-2
- Getting Started..... E-3

Operation

- Display Symbols..... E-4
- Replacing the Batteries..... E-4

References

- Other Information E-5
- Before Requesting Service E-7
- Maintaining the Scale..... E-7
- Specifications..... E-8

- Instructions in French..... F-1

- Instructions in Spanish..... S-1

Precautions

We strongly recommend that you read this instruction manual carefully prior to using the scale the first time.

Installation

- The scale is designed for indoor use.
- Keep the temperature and the humidity of the installation site within the specified operating temperature range / humidity.
- Place the scale on a solid and level floor.
- Do not place the scale on a soft, wet or slippery floor.

Precautions for use

- Do not use the scale if its glass is cracked or broken.
- Do not step onto your scale with wet feet. Make certain that feet are dry to prevent slipping or serious injury.
- Gently step onto the center portion of your scale. Do not jump on the scale.
- Do not disassemble or modify the scale. Prevent foreign matter (staples, clips, water) from invading the scale.
- Do not place the scale in areas where it can be tripped over or kicked.
- Place the scale where it is not exposed to direct sunlight.
- Place the scale where it is not exposed to external noises and strong electromagnetic waves.
- Place the scale where the radio signal between the scale and ActiLink™ won't be significantly impeded by metal or walls.
- Do not place the scale where there are corrosive, flammable, or explosive substances.
- The maximum weighing capacity of the scale is 330 lbs (150 kg). Do not place anything on the scale that is beyond the weighing capacity.
- Do not attempt to service or repair the scale.

Welcome


Congratulations on purchasing a state-of-the-art A&D Medical UC-324THW/UC-324TMW/UC-324THX Precision Scale — one of the most technologically advanced weight scales available today! As part of the Wellness Connected family of products, the UC-324THW works seamlessly with A&D Medical's other Wellness Connected™ products. Designed for ease of use and accuracy, the UC-324THW will effortlessly integrate into your daily wellness regimen.

Physicians agree that daily self-monitoring is an important step individuals can take to maintain overall health. The Wellness Connected family provides a complete system for conveniently tracking vital health information. Use a Wireless Activity Monitor with your new Wireless Scale, or pair it with a Wireless Blood Pressure Monitor, or both!

A&D Medical has been manufacturing high quality healthcare and wellness products for over 20 years. We are committed to providing you and your family with monitoring devices that provide the utmost in accuracy and convenience.

A&D Medical – Designed for Life

Scale Components

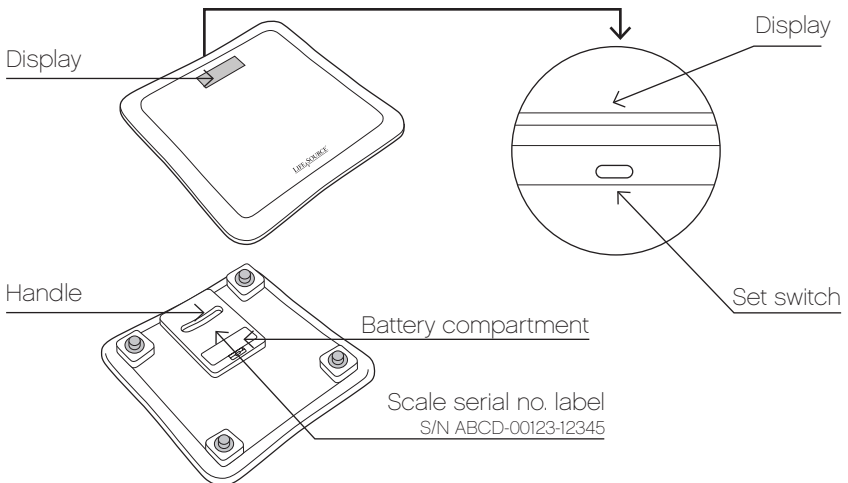Getting Started

Your UC-324THW Wireless Precision Scale will work as a stand-alone device. To use with Wellness Connected follow the below steps.

Before plugging the ActiLink™ USB Transceiver into your computer:


1. Go to www.wellnessconnected.com/getstarted
2. Choose the service plan that is best for your needs.
3. Purchase that plan.
4. Download the Wellness Connected software that is appropriate for your computer.
5. Run the software and follow the installation process.
6. Enter in the required fields (your health information and product serial numbers).
7. Choose your widgets.
8. Set your goals.
9. Start using Wellness Connected!

Using your UC-324THW Precision Scale

1. Insert batteries and place your scale on a solid, level floor. You may change between pounds and kilograms with the Set Switch.
2. Tap the scale with your foot to activate. Wait until the display reads 0.0 lb (0.0 kg).
3. Step on the scale gently and stand still during measurement. Your weight is displayed after the O symbol appears.
4. Step off the scale. The scale automatically transmits your weight wirelessly if your scale is within range of an active ActiLink USB transceiver. If the scale is not in range, the data will be stored in memory until the scale is activated in range of an active ActiLink. The UC-324THW can store up to 63 measurements.

Important note – The ActiLink USB Transceiver will not work when plugged into a computer in sleep, suspend, or hibernation mode.

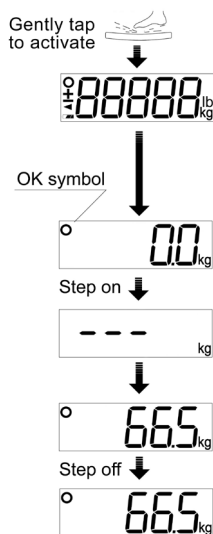Important note – The range of the UC-324THW may be affected by walls, materials around or between the UC-324THW and ActiLink, and other environmental factors.

Display Symbols

Symbols	Description	Operation
	This symbol appears when the measurement display is stable.	The measurement is complete, you may step off.
	The progress bar displays when the scale is finding a stable measurement. Up to five dashes will appear before a stable measurement is reached.	Continue standing still on the scale until a stable measurement is displayed.
	The battery capacity is very low.	Insert new batteries.


Replacing the Batteries

Caution

- Insert the battery with the positive (+) and negative (-) terminals properly aligned with the corresponding symbols in the battery compartment.
- Replace all four batteries with new ones, when  symbol appears.
- Do not mix used and new batteries. It may result in damage to the scale or shortened battery life.
- When the batteries are depleted, nothing is displayed even if the scale is tapped to activate.
- Remove the batteries from the battery compartment when the scale is not to be used for a long time. They may leak and damage the scale.

Steps for Replacing Batteries

- 1 Remove the battery compartment cover on the underside of the scale by gently pushing the hook in the direction indicated by the arrow and lifting the cover backward.
- 2 Install four new batteries with the positive (+) and negative (-) terminals properly aligned with the corresponding symbols in the battery compartment.
- 3 Replace the cover by sliding it into the compartment and gently pressing into place.

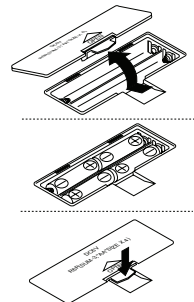Other Information

Our products are designed and manufactured using the latest scientific and technological methods, and offer accurate and easy-to-use home monitoring and treatment options. In conjunction with our UC-324THW Wireless Precision Scale, don't forget to try our other Wellness Connected wireless devices including the Wireless Automatic Blood Pressure Monitor (UA-851THW) or the Wireless Activity Monitor (XL-20). Our full product line includes: Digital Blood Pressure Monitors, Manual Blood Pressure Monitors, Blood Pressure Cuffs, Stethoscopes, Personal Health Scales, and Digital Thermometers.

Warranty:

This Wireless Precision Scale is covered by a **2 year** warranty.

A&D Medical (the "Company") warrants to the first retail purchaser that this A&D Medical Wireless Precision Scale (the "Product") will be free from defects in material and workmanship for two years from the time of purchase. To obtain warranty service for such defects, return the Product, freight and insurance prepaid, together with satisfactory proof of date of purchase such as a sales receipt and shipping costs to:

A&D Medical
Attn: Warranty Department
1756 Automation Parkway
San Jose, CA 95131

Include a check for return shipping and insurance. Contact A&D Medical at 1-888-726-9966 for appropriate shipping costs.

At A&D Medical's option, the Company will repair or replace any Product found to be defective and return it to you. (Any replacement product may be new or refurbished.) This warranty will not apply if the defect or malfunction results from: subjecting the Product to any but the specified voltage; modification, alteration or repairs of the Product by persons not authorized by A&D Medical; misuse, abuse, damage during transit or lack of reasonable care with respect to the Product; or failure to follow the written instructions enclosed with the Product.

THE WARRANTY AND REMEDY PROVIDED ABOVE ARE EXCLUSIVE AND IN LIEU OF ALL OTHER EXPRESS WARRANTIES. IN NO EVENT SHALL A&D MEDICAL BE LIABLE FOR ANY LOSS, INCIDENTAL, OR CONSEQUENTIAL DAMAGES, RESULTING FROM THE USE OR INABILITY TO USE THE PRODUCT, WHETHER RESULTING FROM BREACH OF WARRANTY OR ANY OTHER LEGAL THEORY, EVEN IF A&D MEDICAL HAS BEEN ADVISED OF THE POSSIBILITY OF SUCH INCONVENIENCE, DAMAGE OR LOSS.

Some states do not allow the exclusion of incidental or consequential damage, so that the above exclusions may not apply to you. This warranty gives you specific legal rights, and you may also have other rights that may vary from state to state.

A&D Medical assumes no responsibility with regard to the use of this product. Every effort has been made to ensure that the information in this manual is accurate. A&D Medical is not responsible for printing or clerical errors.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference and (2) this device must accept any interference received, including interference that may cause undesired operations. Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate this device.

Other Information

Disclaimer:

A&D Medical assumes no responsibility with regard to the use of these products. Every effort has been made to ensure that the information in this manual is accurate. A&D Medical is not responsible for printing or clerical errors.

These devices comply with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference and (2) this device must accept any interference received, including interference that may cause undesired operations. Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the devices.

Contact Information:

For more information regarding use, care, or servicing of your UC-324THW Precision Scale, contact:

A&D Medical

A division of A&D Engineering, Inc.

1756 Automation Parkway

San Jose, CA 95131

A&D Medical (Toll-Free): 1-888-726-9966

www.andmedical.com

Before Requesting Service

Problem	Recommended Action
Nothing appears on the display, even when knocking the scale or pressing the switches.	<ul style="list-style-type: none">• Reinstall the batteries with negative and positive terminals matching those indicated on the battery compartment.• Replace all batteries with new ones.
<i>E-0</i> or <i>E-</i> s displayed.	<ul style="list-style-type: none">• Is anything touching the scale?• Is there anything under the scale?• Remove and reinstall the batteries. Use new ones.
<i>888.8</i> is displayed and the scale is turned off.	
There is a large difference between actual weight and measured weight	<ul style="list-style-type: none">• Place the scale on a solid, firm, and level surface.• Did you allow the scale to zero before applying weight?• Is there anything under the scale?
ActiLink/computer software is not receiving new measurements	<ul style="list-style-type: none">• Is your computer on and your ActiLink plugged in?• Move the scale closer to the ActiLink with a clear line-of-sight and take a measurement.• Replace all batteries with new ones.
<i>H</i> is displayed	The 330lb (150kg) weight capacity has been exceeded. Please reduce weight and re-weigh.
<i>L</i> is displayed	Re-weigh and wait until weight is displayed on the screen and/or a final beep is heard.

Note: If the actions described above do not solve the problem, please contact A&D Medical. Do not attempt to open or repair this product, as any attempt to do so will make your warranty invalid.

Maintaining the Scale

- Clean the scale with a soft dry cloth.
- Never use water, detergent, alcohol, benzene, thinner or other harsh chemicals to clean the scale.
- Remove the batteries if the scale is not going to be used for a long time.

Specifications

Measuring method	Electric resistance
Digital display	Liquid crystal display
Maximum capacity (MAX)	330 lbs / 150 kg
Minimum display (d)	0.2 lbs / 0.1 kg
Power source	Four type AA batteries (1.5 V)
Dimensions	13.8 (D) x 13.8 (W) x 1.5 (H) in 350 (D) x 350 (W) x 39 (H) mm
Weight	Approximately 4.8 lbs (2.2 kg), including the batteries
Accessories	Instruction manual, batteries (AA x 4)

Note: Specifications are subject to change without prior notice or any obligation on the part of the manufacturer.

Important!

If you need assistance with the set-up or operation, we can help! Please call us **FIRST** before contacting your retailer at A&D Medical Customer Service Line 1-888-726-9966 (Toll Free). A specially trained representative will assist you.

A&D Medical
A division of A&D Engineering, Inc.
1756 Automation Parkway
San Jose, CA 95131
www.andmedical.com

© 2010 A&D Medical. All rights reserved.

LifeSource, the LifeSource logo, Wellness Connected, and the Wellness Connected logo are trademarks of A&D Medical registered in the U.S.

The UC-324THW Precision Scale and ActiLink USB Transceiver are covered by one or more of U.S. patent numbers 6336365, 6018705, 6052654, 6493652, 6298314, 6560903, 6611789, 6536139, 6876947, 6882955, 7187924 and patents pending.

SPECIFICATIONS SUBJECT TO CHANGE WITHOUT NOTICE.

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Una división de A&D Engineering, Inc.