
**Digital Thermometer
with Lighted Display**

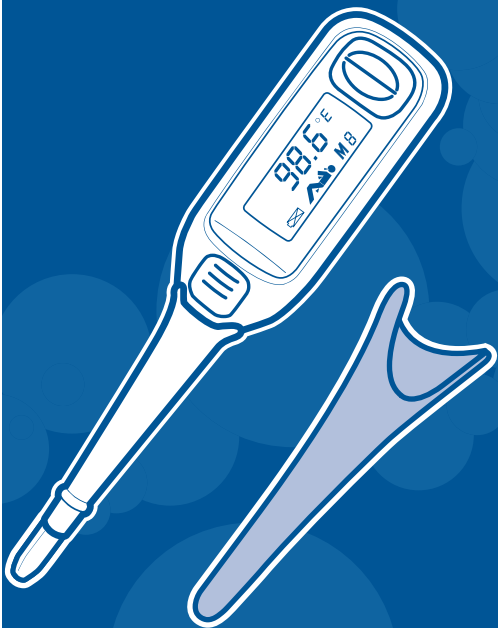GRACO®

OWNER'S MANUAL

KEEP INSTRUCTIONS FOR FUTURE USE.

©2010 Graco PD135573A 3/10

WARNING

Failure to follow these warnings and the assembly instructions could result in serious injury or death.

- **WARNING:** This product is not a toy.
- **DO NOT** immerse any part of the Thermometer in water.
- **DO NOT** place unit near water and moisture. Do not use the Thermometer outdoors. Do not use the Thermometer near possible wet areas, such as a bathtub, shower, wash basin, sink, laundry tub, swimming pool, wet basement, etc.
- **Keep Thermometer away from heat sources**, such as fireplaces, radiators, stoves, and TVs or other appliances. Heat can damage the case or electrical parts.

- **DO NOT** use the Thermometer if:
 - It has been exposed to liquids.
 - It has been dropped or is damaged.

- **DO NOT OPEN THE THERMOMETER.**

No user-serviceable parts inside.

- Do not drop the thermometer!
Protect it from severe impact and shock.

FOR SAFE BATTERY USE:

Keep the battery out of children's reach. Any battery may leak battery acid if mixed with a different battery type, if inserted incorrectly (put in backwards) or if battery is not replaced or recharged at the same time.

Any battery may leak battery acid or explode if disposed of in fire or an attempt is made to charge a battery not intended to be recharged. Never recharge a cell of one type in a recharger made for a cell of another type.

Discard leaky battery immediately. Leaky battery can cause skin burns or other personal injury. When discarding battery, be sure to dispose of them in the proper manner, according to your state and local regulations.

Remove battery before storing product for a prolonged period of time. Battery left in the unit may leak and cause damage.

Recommended type battery—
1 disposable, size CR1225 3-volt lithium battery. Change battery when product ceases to operate satisfactorily.

⚠️ WARNING

Modifications not authorized by the manufacturer may void user's authority to operate this device.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules.

These limits are designed to provide reasonable protection against harmful interference in a residential installation.

This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:


- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

- (1) This device may not cause harmful interference, and
- (2) This device must accept any interference received, including interference that may cause undesired operation.

INTRODUCTION

IMPORTANT: The use of this thermometer should not replace a visit to your doctor. It is recommended that you advise your doctor any time a temperature is above 99.5°F (37.5°C.)

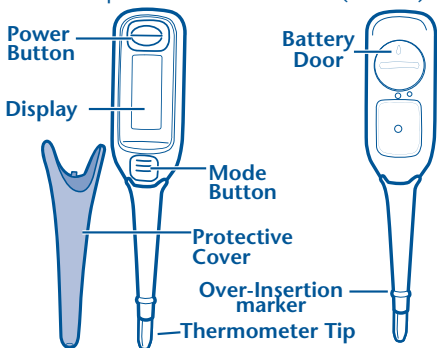ORAL, RECTAL, AND AXILLARY (UNDERARM) TEMPERATURE RANGES

Normal body temperature varies from person to person and depends on factors such as age, gender, physical activity, environmental temperature, and time of day. Normal oral temperatures may range from 95.9°F (35.5°C) to 99.5°F (37.5°C). An oral temperature can also be as much as 0.9°F (0.5°C) lower than a rectal reading and 0.9°F (0.5°C) higher than an axillary reading.

When using after a bath:

- If reading is lower than 89.6°F (32.0°C) the LCD will display 4 dashed lines **A**.
- If reading is higher than 109.3°F (42.9°C) the LCD will display an "H" **B**.


- If temperature is higher than 99.5°C (37.5°C), you will hear 10 beeps and the LCD will flash, letting the user know that a fever has been detected.
- Wait 30 minutes after physical exercise, bathing, eating, or drinking before taking a temperature reading.
- Do not expose the thermometer to extreme temperatures, direct sunlight, dust, or electricity.

BATTERY INSTALLATION

The thermometer includes one CR1225 3-volt lithium battery already installed. The approximate life of this battery is 3 years. When the low battery symbol (⊠) appears on the LCD, follow these instructions for changing the battery.

1. Use a coin to open battery cover from the back of the thermometer ①-②.
2. Remove the battery using a non-metallic tool and discard immediately ③.
3. Install a new CR1225 3-volt lithium battery with the positive side (+) facing up ④.
4. Replace the battery cover ⑤.

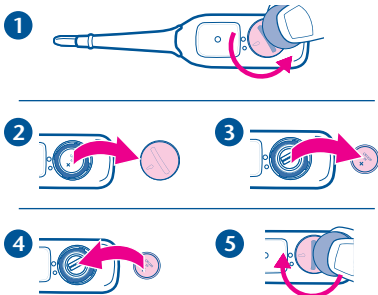
IMPORTANT: Only use a CR1225 3-volt lithium battery. Use of another type of battery may present a risk of fire. The battery may explode if it is handled improperly. Do not recharge, disassemble, heat the battery above 212°F (100°C), or incinerate. Dispose of all used batteries in the appropriate way.

INSTRUCTIONS

For First Time Use

Please follow these instructions to insure the most accurate reading possible.

Remove and discard the protective cover on the front of the thermometer.

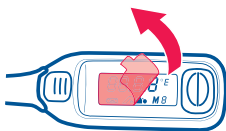How to use:


1. The thermometer default measurement scale is Fahrenheit (F). To change the scale to Celsius (C), press the power button and hold until the current scale setting appears. Release and press the power button again to change the measurement scale. After 3 seconds the thermometer will enter the "ready for measuring" mode automatically.

To select your preferred mode:

2. Turn off the thermometer. Press the Mode button and release. Cycle through the 4 modes by pressing the Mode button until you find your preferred mode.

Each mode will point to a different body part and the stand by icon will display the full body as shown below. After 8 seconds the thermometer will enter the “ready for measuring mode automatically, or you can turn off the thermometer and turn it back on again and you will hear a beep, the thermometer will flash all of the display icons.

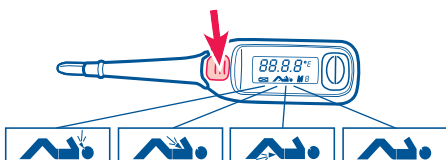NOTE: After the first use, the prior reading and a figure indicating the last measurement method will automatically appear for 2 seconds. Each time you turn the thermometer on, it will automatically default to the same method as the prior reading.

For Oral Measurement


The thermometer is ready to use when the body part symbol is flashing.

Keep your mouth closed for 5 minutes before taking the temperature. Place the tip under the tongue at the back of the mouth. For the most accurate read, place the tip as far back under the tongue as possible. Press the tongue down over the tip, and close your mouth. Hold the thermometer steady keeping the tip in the same position for 5 seconds. After approximately 5 seconds you will hear a five short beeps indicating the measurement is complete.

- (3 X “bi-bi-bi”) Temperature is higher than 42.9°C.
- (10 X “bi-bi-bi-bi”) Fever is detected.

NOTE: If the tip is not placed in the proper position, the temperature reading can be as much as 1.2°F (0.7°C) lower than one taken in the correct position.


- ✓ = correct areas
- X = incorrect areas

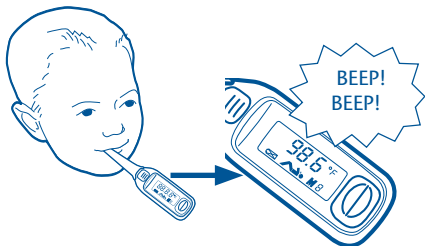
For Axillary (Underarm) Measurement


The thermometer is ready to use when the body part symbol is flashing.

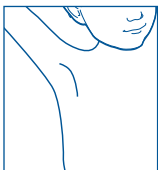Press the Mode button until you see the Axillary measurement icon. Once the icon is showing, allow the thermometer to stay idle for 8 seconds before taking a measurement.

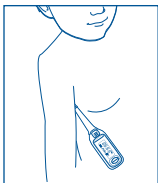Make sure the underarm is dry and there is no material between the side and the arm. Hold the arm close to the side for 5 minutes before taking a temperature reading. Point the thermometer upward and place the tip well in the underarm.

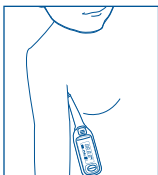Hold the thermometer with the tip slightly bent and the arm firmly pressed against the side for approximately 10 seconds. You will hear a 5 beeps when the temperature reading is complete.


For Rectal Measurement

The thermometer is ready to use when the body part symbol is flashing.

Lubricate the tip of the thermometer with a water-soluble jelly for easier insertion.

Press the Mode button until you see the rectal measurement icon. Once the icon is showing, allow the thermometer to stay idle for 8 seconds before taking a measurement.

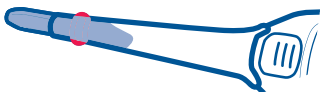Gently slide the tip no more than 3/4 inch (19 mm) into the rectum. DO NOT push thermometer past the over-insertion marker.

Hold the thermometer in place for 5 seconds. After approximately 5 seconds you will hear five short beeps, indicating the measurement is complete. If you want to use the thermometer in another mode after using it rectally, be sure to clean it very thoroughly with isopropyl alcohol.

NOTE: If you detect any resistance when inserting the thermometer, stop immediately.

Stand by Icon

This is a non use mode. When you reach this symbol, you have cycled through all the modes available. Press the Mode button again until you have reached your preferred mode.

3. Turn off the thermometer after use by pressing the Power button. If you forget to turn it off it will automatically shut off after 3 minutes.

4. Clean the tip with isopropyl alcohol before and after each use.

5. Replace the protective cover after each use and keep it on when storing the thermometer.

NOTE:

- Do not immerse the thermometer in water.


MEMORY FUNCTION

Using the Memory function

With the thermometer turned off, press and hold the Mode button for at least 5 seconds until the "M" icon appears on the LCD display.

Pressing the Mode button again will provide you with the last temperature reading followed by a beep. If you continue to press on the Mode button, you can cycle through, up to 10 readings. Each reading will display a number from M0 to M9.

NOTE: To turn off, press the Power button or leave thermometer idle for about 3 minutes.

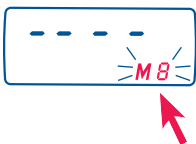CARE AND CLEANING

1. The thermometer is a delicate electronic device. Although the probe can be cleaned with isopropyl alcohol, the body of the thermometer should never be submerged in water or any other liquid.

2. Clean the probe with a cotton swab moistened with isopropyl alcohol after every measurement and before each use. Clean the body of the thermometer with a soft dry cloth. Never use an abrasive.

3. Avoid touching the probe with your fingers.

4. Do not expose the thermometer to direct sunlight.

5. Do not disassemble the thermometer.

6. Use only a CR1225 3-volt lithium battery.

This thermometer is not a toy. Do not allow children to play with it. Keep out of reach of children when not in use.

TROUBLESHOOTING

Display	Possible Cause	Solution
	When using after a bath: Measured temperature is too high 109.3°F (42.9°C),	Take a new reading and wait for the beep before removing the thermometer.
	When using after a bath: If reading is lower than 89.6°F (32.0°C)	Take a new reading and wait for the beep before removing the thermometer.
	The thermometer needs to be reset or the battery should be replaced.	Reset the thermometer by removing and reinstalling the same battery; replace the battery with a new one; call customer service.
	Unstable temperature reading.	Take a new measurement while holding the thermometer in a stable position.
	When the low battery symbol (☒) appears on the LCD,	Replace old battery with a new one immediately.

TECHNICAL SPECIFICATIONS

- For home use only
- Temperature readings are to the nearest 0.1°F/C
- Temperature measurement range: 89.6°F-109.3°F (32.0°C-42.9°C)
- Operating temperature range: 50°F-104°F (10°C-40°C)
- Storage Environment: -25 °C to 60°C, humidity of 15 % to 95 % (noncondensing)
- Requires one CR1225 3-volt lithium battery (included)

Applied Standards and Approvals

Device standard: Device corresponds to the requirements of the standard for electronic thermometers ASTM E 1112 - 98 (Accuracy).

Electromagnetic Compatibility: Device fulfills the stipulations of the standard IEC 60601-1-2.

Replacement Parts Warranty Information (USA)

To purchase parts or accessories or for warranty information in the United States, please contact us at the following:

Graco Children's Products
3 Glenlake Parkway
Atlanta, GA 30328

www.gracobaby.com
or
1-800-345-4109

Product Registration (USA)

To register your Graco product from within the U.S.A. visit us online at

www.gracobaby.com/productregistration.

We currently do not accept product registrations from those living outside the United States of America.

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