

DNS REVO (DNS80)

User Manual



EACH DNS REVO WITH ONE LED LIGHT



EACH PACKAGE MATCH TWO ROLLERS

Instruction


Micro-needle therapy or Derma roller is used for skin micro roller "minimally invasive" stimulation, in a short time to create a large number of micro-channels in the epidermis and dermis, subcutaneous tissue formed between a superior delivery system, the activity of the product ingredients and nutrients to be quickly absorbed through the skin, greatly enhance the product's efficacy.

"DNS REVO" more with the ***optical micro-needle therapy, vibrating massage therapy, bionic and galvanic energy therapy three power technology***, compared with the ordinary micro-needle, "DNS REVO" has a quick, less invasive, healing time is short, and according to skin surface, middle, or deep skin and hypodermic layer, then processing selective care; and realize painless treatments.

Parts Instruction



1. Current Intensity Display Light
2. Vibration Function Start/Pause
3. Galvanic Intensity Press Button
4. Galvanic Weaken Press Button
5. Roller Protecting Cover

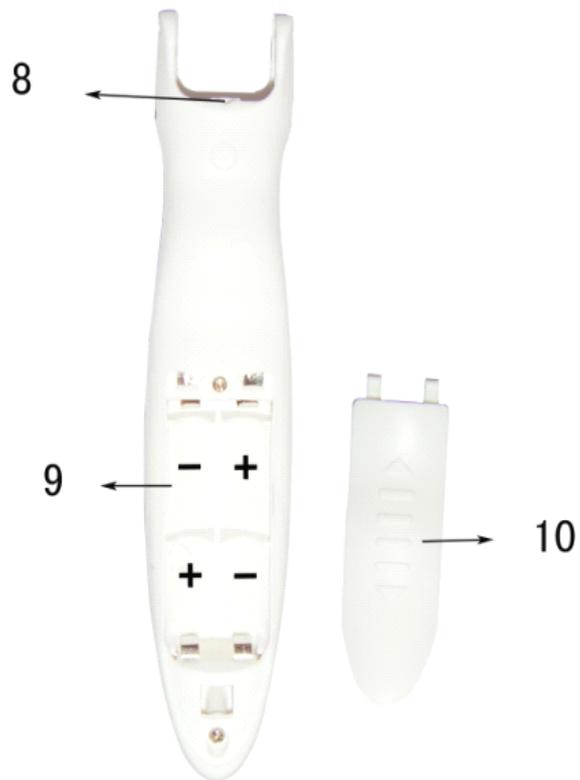

6. DNS BIO Galvanic Roller Head
7. Roller mounting board
8. LED therapy
9. Double AAA Battery Position
10. Battery Cover


REMARKS

DNS REVO without match the battery, yours can buy it your Local Place, it is very easy to handle it. In this way, it fits international Environmental Protection Standards, also easy to take and use etc.

Basic Operation of Instrument

1. Insert two AAA batteries



The "positive and negative" pole should accordance to the instruction of handle.

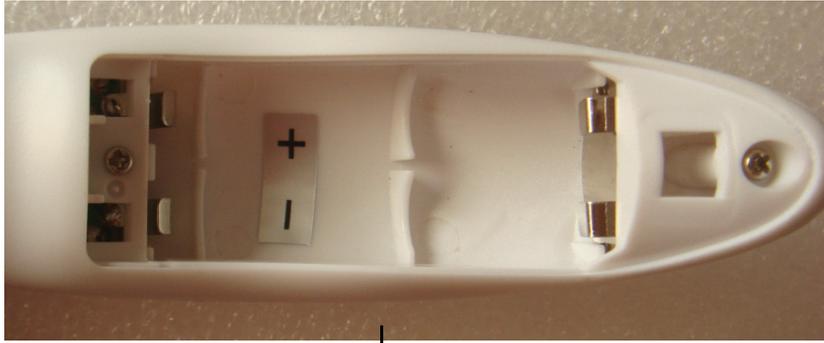

Put the cover on the pathway, slide in slowly

2. Select the desired roller head and assemble it into the handle, then remove the cover.



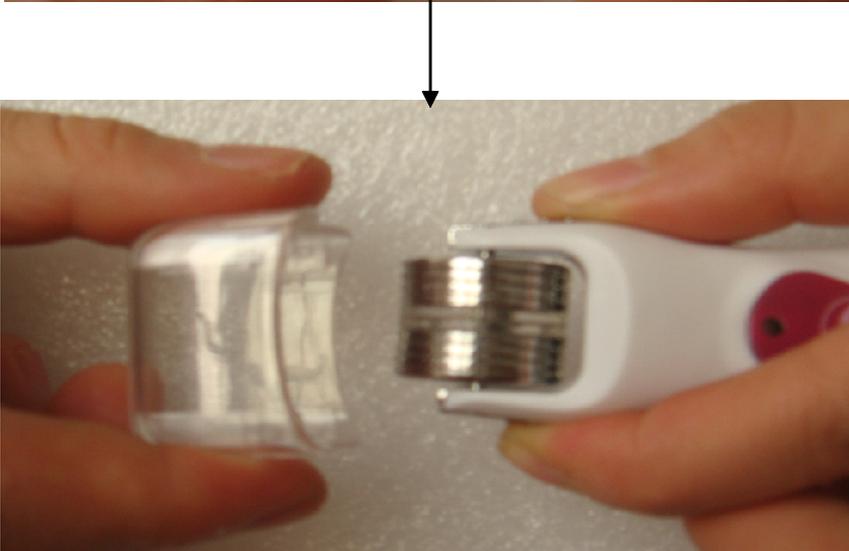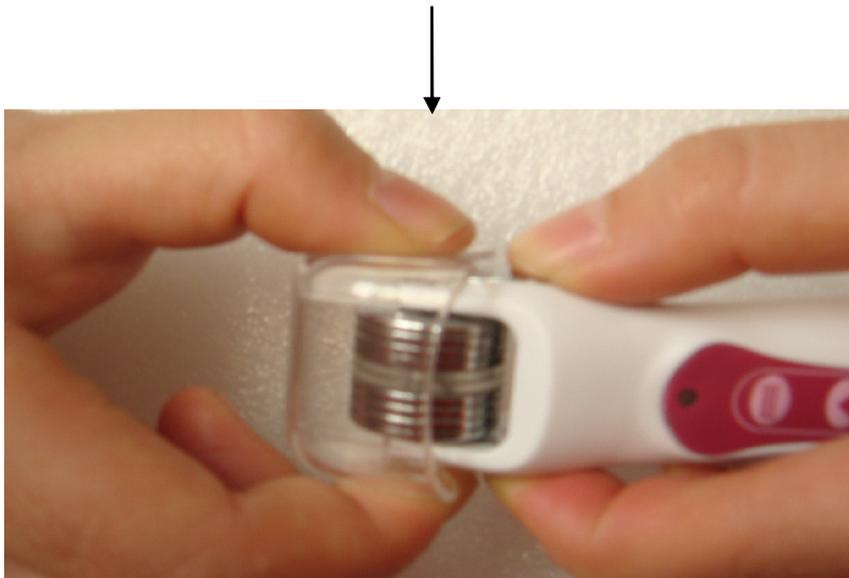

3. Start the optical led therapy

Press "3" one time, the "1" turn green, at the same time, the "8" led light turn on, the LED therapy start now



4. Adjust the electric energy bionic intensity.

After the last step. The micro-current function turned on together with the led therapy, but it's in a low intensity, there are total 6 channel for the intensity of micro-current,

Channel 0: the original state

Channel 1: press 1 times, the "1" indicator turn light (dark green)

Channel 2: press 2 times, the "1" turn more lighter (bright green)

Channel 3: press 3 times, the "1" turn red, (dark red)

Channel 4: press 4 times, the "1" turn lighter (bright red)

Channel 5: press 5 times, the "1" shows little difference with channel 4, but the intensity will be the strongest

(To turn off the micro-current, press 4 to channel 0)



5. Starts vibrating massage function.

Press "2" to start skin massage function, there are three channel for this function.

Channel 0: the original state

Channel 1: it's more appropriate facial care

Channel 2: arm or legs care


Channel 3: for stomach or around

Note: the intensity of vibrating becomes strong from channel 1 to channel 3

When start the vibration function, the led therapy start the same time



6. Treatment on Skin:

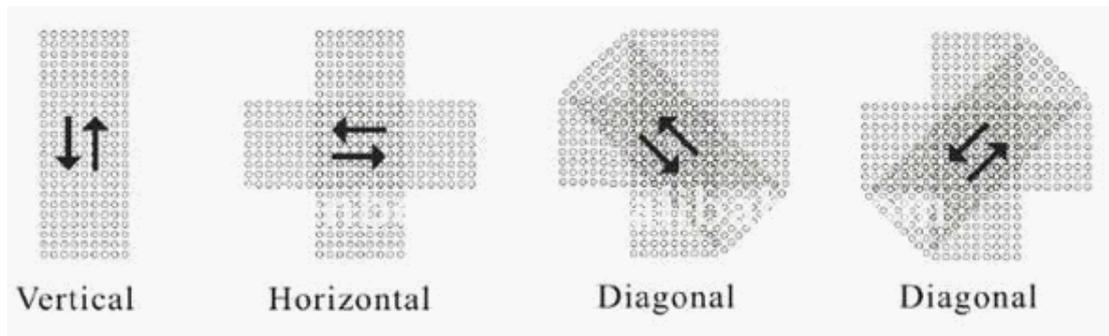

Note: only for your reference

Principle and Effectiveness

First, Optical Therapy

The theoretical basis of Optical Therapy is biological tissues absorb light energy and light energy into heat energy and chemical energy, resulting in the human body produce a series of chain reactions, summed up in four types: light-induced decomposition of the light-induced oxidation, light-induced sensitization of polymerization and light, produced by different wavelengths of light have different biological effects:

1. 405nm blue light to eliminate acne inflammation, anti-bacterial anti-inflammatory, inhibiting the secretion of sebum, anti-aging effects, but also make your skin soft and smooth, white charming!

Scope: Treatment of acne, acne, acne and other skin.

2. 633nm high-energy narrow-band red light can accelerate the healing of minimally invasive surface; the human muscle and skin cells to 5 times than the normal rate of growth, increased endurance and resistance to fatigue, promote cell metabolism.

Scope: anti-aging, wrinkle repair and restore skin elasticity.

3. 590nm yellow supplements cellular energy, and promote the role of the glands, digestive aid, treatment of skin diseases, strengthening the immune system.

Scope: shrink pores, dull skin repair, improve rough skin, redness and other skin.

4. 518nm with a balance of green cured swollen detoxification, calm the skin

effect.

Scope: oily skin, spots and other pigment felt like skin.

Second, vibration massage function:

Mechanical vibrations acting on the skin will change to physical vibration, vibration massage can promote blood circulation and lymphatic circulation, promote skin metabolism and activation of cell movement, vibration massage can eliminate skin fatigue, relax tense skin, so that micro-needle roller acting on the skin more uniform intensity.

Third, the bio-electric energy:

Bio-electric energy, which is very close to the body's biological micro-current, and its main role is to help restore damaged and aging cells to normal, to maintain normal cell metabolism and circulation, improve the ability to repair damaged cells and promote aging of cells gradually Reply to a young man of energy, activate cells, dormant cells active again; accelerated collagen and elastic fibers and collagen fibers, healing ability, to achieve a deep nourish and improve the results.

Specification

One. 0.2MM, 0.25MM, 0.3MM fit for the skin too thin, sensitive skin, or redness of skin collagen activate the repair, maintenance and daily skin care routine to help the absorption of nutrients.

Scope: whitening, moisturizing, light spot etc.

Two. 0.3MM and 0.5mm fits Stratum epidermis which can create many microtubules of penetrate cells, then form a high-quality nutrient delivery systems, micro-needle can penetrate to the dermal layer of nutrients, can relieve eye fine lines around the scope of such a thin epidermis, elimination of smallpox in India , sunburn and other effects lasting and reliable.

Scope: skin whitening, Dodge stain, improve eye wrinkles, dark circles, bags under the eyes, facial skin tightening and improved.

Three.1.0mm and 1.5mm fits Stratum epidermis which can create many microtubules of penetrate cells, then form a high-quality nutrient delivery systems, micro-needle can penetrate to the dermal layer of nutrients, nutrient uptake capacity of activated fibroblasts. Scope can be whitening, water, enhance facial tissue tightening, treatment of skin, hair loss, acne, hair is thick, scars, uneven and so on.

Scope: for the treatment of deeper wrinkles, cellulite, stretch marks, weight loss, weight, eliminates tattoos, scars and so on.

Four.1.5mm and 2.0mm

Apply massage to penetrate to the dermal layer, by stimulating the skin's natural healing ability to induce proliferation of natural collagen, a significant effect for the loss of wrinkles, stretch marks and also diminish general scarring, bean pit, scope: reduce the appearance of wrinkles, induce bone collagen, treatment cellulite, stretch marks, such as obesity pattern thickness skin problems.

DNS COMPARISON WITH OTHER DERMAROLLERS

FIG ONE



FIG TWO

	192 Pins	540 Pins	DNS - REVO
Punctures	Extensive	Mild	Mild
Pain	Extensive	Mild	Mild
Recovery	Slow	Medium	Fast
Result	Excellent	Fair	Astonishing
Life Span	Short	Fair	Long

CLINICAL FEEDBACK



Before Treatment



After Treatment

Notes

1. After care, it is not exposure to high temperature, heat the place, avoid saunas, hot springs for a long time within three days.
2. Please don't use the scrub, exfoliating products after nursing; please do not make diamond dermabrasion or micro-dermabrasion treatment.
3. To maintain general skin care, do not use irritating skin care products.
4. Please don't operate it on the wound and inflammation of the skin.
5. It should from weak to stronger when you adjust the current or galvanic intensity.
6. Please don't use it on the eye lid, Adam's apple and also the part of heart.
7. Each roller allows only one guest use, and also must be sterilized before use.

Forbidden of people



1. infectious, acute illness
2. heart disease, cardiac pacemakers and devices people
3. severe hypertension, cancer disease, asthma
4. bleeding disorders or trauma, dermatitis, skin
5. the menstrual period not to use care in the abdomen
6. pregnant women
7. the medical use of a silicone plastic parts
8. implanted metal