
**Safety 1<sup>st</sup>**

# Accu-Touch Forehead Thermometer

Model 49529
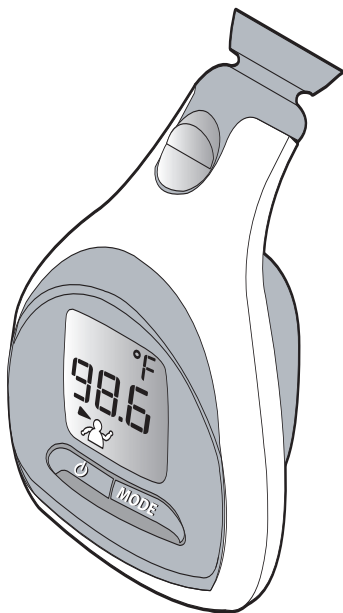
## User Guide

---

Questions? Comments?  
Missing Parts? Contact Us:

[www.safety1st.com](http://www.safety1st.com)  
(800) 544-1108

©2004 Dorel Juvenile Group  
Made in ISRAEL  
[www.safety1st.com](http://www.safety1st.com)  
(800) 544-1108

Styles and colors may vary.  
Distributed by Dorel Juvenile Group  
2525 State Street,  
Columbus, IN 47201-7494  
Dorel Sales and Distribution Canada  
873 Hodge St., St-Laurent,  
QC Canada H4N 2B1  
06/05 4358-3760

# WARNING

**BEFORE USING THIS THERMOMETER, PLEASE READ CAREFULLY:**

- This thermometer is not a toy and should be kept out of a child's reach when not in use.
- Keep thermometer away from hot or boiling water, chemicals, and direct sunlight.
- Avoid sudden shock and dust.
- Make sure the thermometer has been at room temperature for 15 minutes before use.
- Remove batteries if stored for a long period of time.
- Do not disassemble.

## **SAFETY INFORMATION:**

- High, prolonged fever requires medical attention especially for young children. Please contact your physician.
- Do not allow children to take their temperatures unattended.
- Read these instructions carefully to ensure accurate temperature readings. Temperature readings are affected by many factors including exercise, drinking hot or cold beverages prior to taking a reading as well as technique.
- Do not allow children to walk or run during temperature taking. The probe could present a hazard.

## **Introduction**


---

Thank you for purchasing the Accu-Touch Forehead Thermometer, the thermometer that enables accurate temperature measurements of both forehead and underarm in approximately 7 seconds.\*

\*Timing results may vary.

## How Does it Work?

This thermometer is based on R.A.T.E™ (Rapid, Accurate, Temperature, Establishment) technology. The accurate measurement is obtained by using a unique probe, which measures the heat flow from the blood vessels to the skin and converts it to body temperature.
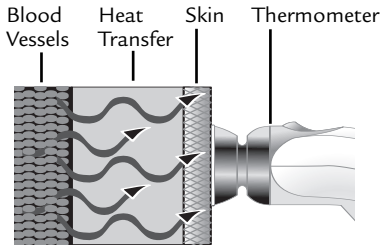
## Body Temperature

What is “normal” body temperature? “Normal” body temperature varies from person to person and changes according to time of day and body location. Hot or cold baths, physical exercise, smoking and drinking hot or cold drinks may change body temperature by one degree or more.

The following table shows the temperature ranges:

	What is fever?	Average normal temperature ranges
Forehead	Over 100.0°F / 37.8°C	97.0°F - 100.0°F / 36.1°C - 37.8°C
Rectal*	Over 100.4°F / 38°C	97.3°F - 100.4°F / 36.3°C - 38°C
Oral *	Over 100.0°F / 37.8°C	97.0°F - 100.0°F / 36.1°C - 37.8°C
Underarm*	Over 99.0°F / 37.2°C	95.4°F - 99.0°F / 35.2°C - 37.2°C

\* Fever definitions according to the American Medical Association

The thermometer is designed to precisely measure the temperature of the forehead between the eye socket and the hairline. It is important to know that the thermal characteristics of each measurement body location (temple, underarm, oral or rectal) are unique. The variance in temperature measurement in different

# Introduction

---

## Body Temperature Continued

body locations is due to the human body and not the thermometer. The Accu-Touch Forehead thermometer readings are equivalent to oral measurements.

### IMPORTANT:

- It is recommended that you measure the body temperature of each individual when healthy to establish a “normal” temperature.
- Due to the unique body physiology of each person, it is recommended that you take a persons temperature at the same temple location to compare measurements.

## Information on Fever

A fever is a body temperature that exceeds an individual’s normal body temperature range. Except in the cases of newborn infants, the presence of a fever is not necessarily a cause for concern. In fact, a fever can help the immune system to work more effectively. However, a physician should be consulted if a fever is accompanied by vomiting; diarrhea; changes in appetite, breathing, or activity level; or, in the case of children, irritability, lethargy, or sleepiness.

## Memory Function


The Accu-Touch Forehead Thermometer memory function stores the last measurement. This measurement will appear for approximately two seconds after turning the thermometer on and the thermometer initializes. The ready display will then show after last reading appears. This feature allows you to compare your most recent measurement to the prior one.

## Automatic Shut Off


Automatic shut off occurs after one minute of no activity.

## Parts Names

**NOTE:** The names of the parts are shown below for your convenience.
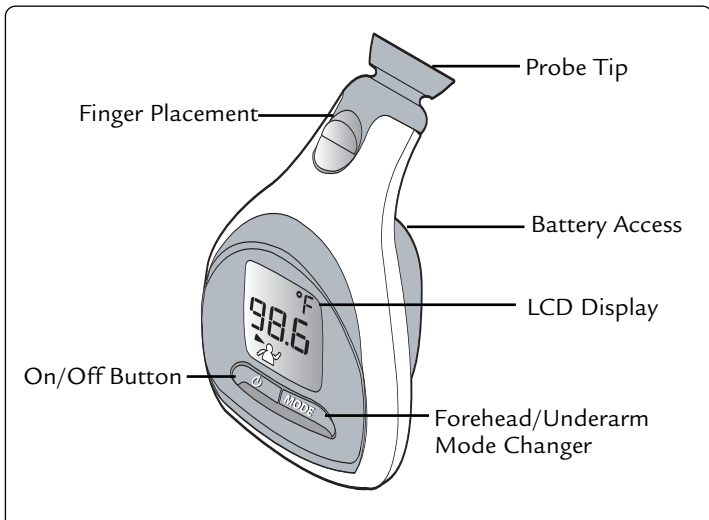
## Thermometer LCD Display

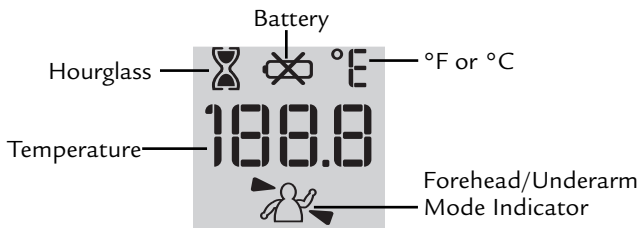**NOTE:** Display is shown with all segments lit.

# Using Your Thermometer

---

**Before You Begin:** What to know before you get started:


**IMPORTANT:** Follow these guidelines:

- Brush hair away from the temple area.
- Press the Mode button to select forehead or underarm mode. Thermometer will start in last mode used.
- Use this thermometer on a clean, dry forehead or underarm. Wipe away any perspiration prior to use.
- Do not try to take a measurement until the two beeps are heard.
- Make sure that the thermometer and the patient do not move while measurement is in progress.
- Correct placement on the temple and good contact with the body ensure reliable measurements.
- As soon as the thermometer is in place, the measurement starts automatically. A running circle will appear indicating that measurement is progressing successfully.
- Wait at least 2 minutes before taking another measurement.
- Once the thermometer is “ready to measure”, it will shut off automatically after 1 minute if it has not maintained good contact with the patient’s skin.
- Clean the thermometer after each use (see “Care and Cleaning” section).

# Using Your Thermometer


## To Take a Forehead Measurement:

1. Press and release the ON/OFF button. You will hear one beep. All segments will be displayed (Figure 1).

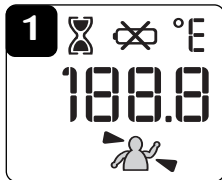
**IMPORTANT:** Press the Mode button to select **forehead measurement mode**. Be sure that an arrow is pointing at the forehead as shown (Figure 2).

2. The result of the previous temperature measurement will be displayed for 2 seconds (Figure 2).
3. Wait for the hourglass icon to disappear, then two beeps will be heard to indicate that the thermometer is “ready to measure”.

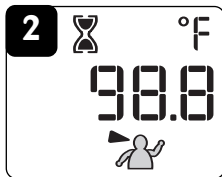4. Wipe away perspiration with a dry towel.
5. Place the thermometer probe in the center of the temple area between the eye socket and the hairline. Gently press it against the temple and ensure good contact with the skin is maintained (Figure 3).

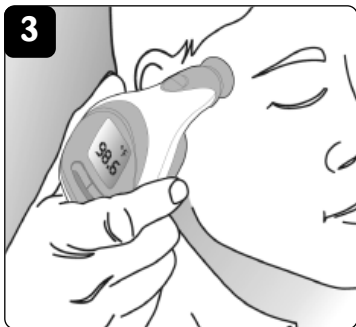6. A temperature reading will be displayed within 7 seconds\*. One long beep will be heard to indicate that measurement is complete and temperature will be displayed for about 1 minute.

7. Press and release the ON/OFF button to turn the thermometer off.

**NOTE:** If you need to take another reading, wait two minutes before taking another measurement. Thermometer will start in last mode used.

\*Timing results may vary.

# Using Your Thermometer

## To Take an Underarm Measurement:

1. Press and release the ON/OFF button. You will hear one beep. All segments will be displayed (Figure 1).

**IMPORTANT:** Press the Mode button to select **underarm measurement mode**. Be sure that an arrow is pointing at the underarm as shown (Figure 2).

2. The result of the previous temperature measurement will be displayed for 2 seconds (Figure 2).
3. Wait for the hourglass icon to disappear, then two beeps will be heard to indicate that the thermometer is “ready to measure”.

4. Make sure the armpit is dry and there is no clothing between the chest and arm. Wipe away perspiration with a dry towel.
5. Place the thermometer probe in the center of the armpit as shown (Figure 3). Gently press the thermometer’s probe in an upward motion.




6. Leave the patient’s arm raised up – **DO NOT BRING THE HAND DOWN DURING THE MEASUREMENT** until completed.

7. A temperature reading will be displayed within 7 seconds\*. One long beep will be heard to indicate that measurement is complete and temperature will be displayed for about 1 minute.

8. Press and release the ON/OFF button to turn the thermometer off.

**NOTE:** If you need to take another reading, wait two minutes before taking another measurement. Thermometer will start in last mode used.

\*Timing results may vary.

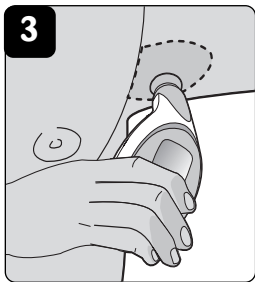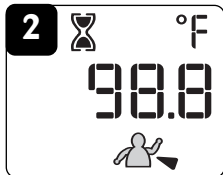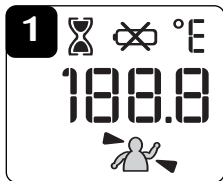
## To Change the Measurement Scale:

**NOTE:** The Accu-Touch Forehead Thermometer can display temperature measurements in either Fahrenheit or Celsius. The unit will be factory preset to Fahrenheit. Use this procedure to change to Celsius or back to Fahrenheit at any time.

1. To switch the display between °F and °C and vice-versa: While off, press and hold the MODE button. Continue holding while you press and release the ON/OFF button (Figure 1).
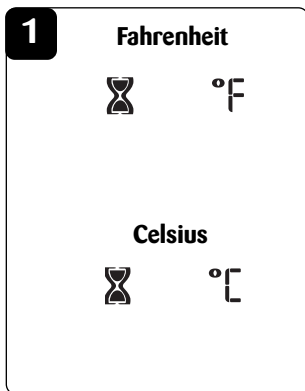
## Care and Cleaning

---

- Wipe the probe with an alcohol swab or cotton swab moistened with 70% isopropyl alcohol before & after each use.
- After cleaning, wait at least two minutes for the unit to adjust to room temperature before taking measurement.
- Use only a soft, dry cloth to clean the body of the thermometer (if necessary).
- Never use an abrasive cleanser or submerge the thermometer in water or other liquids.
- Do not expose the thermometer to extreme temperatures, humidity, direct sunlight or shock. Store at room temperature.


# Battery Replacement

The Accu-Touch Forehead Thermometer includes two AAA batteries already installed which should provide at least 1000 measurements. Approximate battery life is three years. When the Low Battery Display appears, replace batteries before taking a new measurement. 

## Installing and Replacing the Batteries:

The thermometer uses 2 standard alkaline AAA 1.5v batteries. To install/replace :

1. Using a screwdriver, loosen the screw on the battery door and lift battery cover off thermometer.
2. Place batteries in battery compartment, matching polarity according to the diagram.
3. Secure the battery cover using the screw.
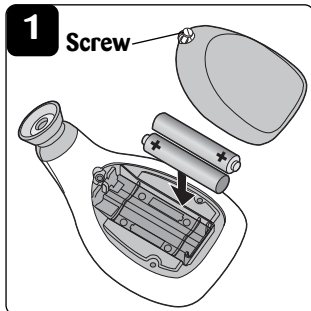
## WARNING

- Remove batteries if stored for a long period of time.
- KEEP BATTERIES OUT OF CHILD'S REACH.  
SWALLOWING A BATTERY COULD BE FATAL.  
BATTERIES SHOULD NOT BE CHARGED OR PLACED INTO EXTREME HEAT AS THEY MAY RUPTURE OR EXPLODE.

## Troubleshooting and Displays

The thermometer has been designed to display feedback/warnings to ensure proper usage and accurate temperature measurement. When an error or indication appears, one “beep” signal will be heard once every 5 seconds for 10 seconds and then the thermometer will turn off automatically. In the event of a malfunction or incorrect temperature measurement, the following warnings will be displayed:

Display	Meaning	Action
	The measured temperature is below 95°F (35°C).	Turn the unit OFF and then ON again. Retake temperature as described in “Using Your Thermometer” section.
	The measured temperature is above 107.6°F (42°C). Or thermometer was placed in body position before two beeps were heard.	Turn the unit OFF and then ON again. To verify the reading, retake temperature as described in “Using Your Thermometer” section.
	The thermometer is too cold or too hot to operate.	Verify room temperature is within 50°F to 104°F (16°C to 40°C). Retake temperature as described in “Using Your Thermometer” section.
	Temperature did not register because of poor contact, patient movement or poor placement.	Turn the unit OFF, wait two minutes, turn it ON again and retake temperature as described in “Using Your Thermometer” section.
	Unit malfunctioned.	Do not use the thermometer - contact distributor for instructions.
	Battery icon is on continuously.	The batteries have enough remaining life for a few more measurements. Replace batteries.
	Only battery icon is on and the device automatically turns itself OFF.	The batteries are drained - replace batteries (see “Replacing Batteries” section).

## Technical Specifications

---

Feature	Performance / Data
Measurement time	7 seconds (timing results may vary)
Measurement locations	Temple artery area of the forehead or underarm
Measurement range	95.5°F to 107.6°F (35°C to 42°C)
Operation environment	60.8°F to 104.0°F (16°C to 40°C) Humidity: 15-95% (non condensing)
Storage environment	-4°F to 122°F (-20° C to 50° C), Humidity: 15-95% (non condensing)
Measurement readout	Fahrenheit or Celsius
Accuracy	Complies with ASTM E 1112-00 standard requirements
Batteries	2 AAA 1.5volt alkaline batteries
Battery life	1000 full duty cycles
Memory	Last reading display when turned on
Auto shut off	After 60 seconds
Probe material	Stainless steel and food grade ABS plastic
Warranty	1 year

# One (1) Year Limited Warranty

---

Dorel Juvenile Group, Inc. warrants to the original purchaser that this product (Accu-Touch Forehead Thermometer), is free from material and workmanship defects when used under normal conditions for a period of one (1) year from the date of purchase. Should the product contain defects in material or workmanship Dorel Juvenile Group, Inc., will repair or replace the product, at our option, free of charge.

Purchaser will be responsible for all costs associated with packaging and shipping the product to Dorel Juvenile Group Consumer Relations Department at the address noted on the front of this document and all other freight or insurance costs associated with the return. Dorel Juvenile Group will bear the cost of shipping the repaired or replaced product to the purchaser. Product should be returned in its original package accompanied by a proof of purchase, either a sales receipt or other proof that the product is within the warranty period. This warranty is void if the owner repairs or modifies the product or the product has been damaged as a result of misuse.

This warranty excludes any liability other than that expressly stated previously, including but not limited to, any incidental or consequential damages.


**SOME STATES DO NOT ALLOW THE EXCLUSION OR LIMITATION OF INCIDENTAL OR CONSEQUENTIAL DAMAGES, SO THE ABOVE LIMITATION OR EXCLUSION MAY NOT APPLY TO YOU.**

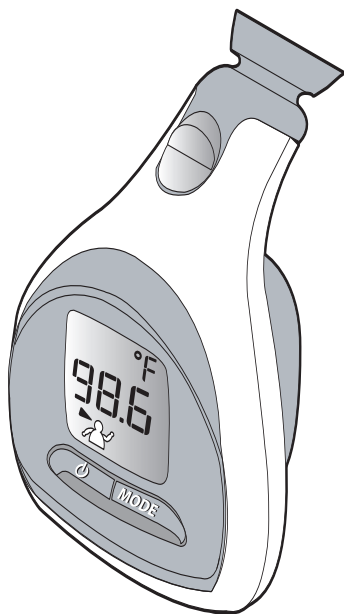**THIS WARRANTY GIVES YOU SPECIFIC LEGAL RIGHTS, AND YOU MAY ALSO HAVE OTHER RIGHTS THAT MAY VARY FROM STATE TO STATE.**

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