




03160

Please read the following instructions and warnings carefully.

Keep this instruction manual for future reference.

For Technical Support Call 1-800-268-6237 or E-Mail at customerservice@summerinfant.com

For use with children from birth and up.

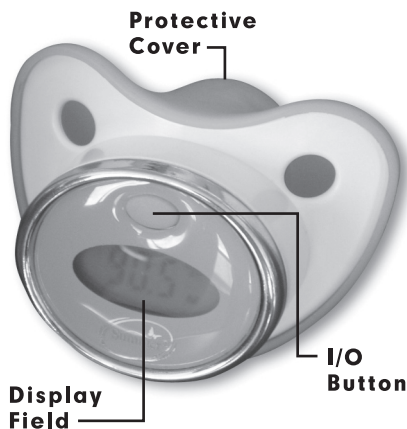

INSTRUCTION MANUAL

PACIFIER THERMOMETER

Congratulations on your purchase of the Pacifier Thermometer. Please read these instructions carefully to ensure accurate temperatures and safe operation. When used properly, this high quality product will give you years of accurate, dependable service.

⚠ WARNING:

DO NOT TIE PACIFIER AROUND CHILD'S NECK AS IT PRESENTS A STRANGULATION DANGER.

IMPORTANT SAFETY INSTRUCTIONS:

NOTE: The Pacifier Thermometer is intended for use with children ages five and under.

1. **THIS UNIT IS A MEDICAL DEVICE** and is not designed for use as a regular pacifier. Store out of reach of children. Do not let children play with the unit. **TEMPERATURE SHOULD BE TAKEN UNDER ADULT SUPERVISION.**
2. Any pacifier wears out. **TO PREVENT POSSIBLE CHOKING HAZARD**, test the strength of the nipple by pulling on the bulb portion before each use.
3. Examine the pacifier bulb before each use for cracks, tears, tackiness or other signs of wear. Discard a worn unit immediately.
4. Components may be harmful if swallowed. Contact your physician immediately should this occur.
5. **FOR ORAL ONLY.** Do not measure temperature at any other site.
6. The thermometer nipple comes in one size (medium) intended to fit most mouth sizes. Please gauge your child's mouth size and the nipple size prior to thermometer use. A child using a nipple that is too large has an increase risk of choking.
7. Do not attempt to disassemble the unit. There are no user serviceable part inside. The BATTERY IS NOT REPLACEABLE. When the battery is almost exhausted, the following symbol will appear in the display window(■□□) indicating the unit should be discarded. Do not use as a regular pacifier after the battery is discharged.
8. Wash the unit in warm soapy water. Do not wash the unit with any thinner or chemical solvent.
9. Do not store unit in direct sunlight or at high temperature.

TEMPERATURE READINGS:

This pacifier is for oral use only. Oral temperature readings are often approximately 1° lower than rectal readings and 1° higher than axillary (under the arm) readings.

Temperature readings vary from child to child. Although there is no singular "normal" temperature reading, a temperature taken orally ranging between 97°F (36°C) and 100°F (37.8°C) is considered normal. Your child's physician should be consulted to determine what temperature range the physician considers to be normal.

The best method to determine your child's normal temperature is to use the thermometer when the child is feeling well. Record readings twice a day (early morning and late afternoon). Take the average of the two temperatures to calculate a normal body temperature.

OPERATING YOUR PACIFIER THERMOMETER:

1. Press the power button. The display will read (188.8^{°F}) indicating that the display is working and a short signaling tone is heard, then the last reading is displayed for 2 seconds, followed by the circuit test. If the circuit is normal it displays (Lo) and a flashing (°F) appear at the upper right of the display. The thermometer is now ready for use.
2. When the temperature reading reaches 89.6°F (32°C) or more, the display will show an actual temperature.
3. Position the pacifier bulb in child's mouth. Note: Drinking or eating prior to use may affect reading. Measurement should be taken one half hour after eating or drinking.
4. Once the degree sign (°) on the display has stopped flashing (approximately one and half minutes), the beeper will sound, meaning the temperature is ready to read. This reading reflects the oral temperature of the child.
5. If the reading is below 99.5°F (37.5°C) the green backlight will be on for 5 seconds. If the reading is equal or above 99.5°F (37.5°C) the red backlight will be on for 5 seconds.
6. It is recommended that the thermometer be left in the child's mouth for 3-5 minutes to allow the nipple to warm up and to obtain a more accurate reading. This is required if the child does not pacify without interruption.
7. The unit will automatically turn off in 3 and half minutes. However, to prolong battery life, press the power button to turn the unit off once the temperature has been taken.
8. Repeat this procedure as necessary if pacifying is interrupted during use.

HOW TO CHANGE TEMPERATURE SCALE:

The Pacifier Thermometer can display temperature measurements in either Fahrenheit or Celsius.

To switch the scale from °C to °F or vise versa, press the O/I button for more than 3 seconds, and the unit is changed from °C to °F, then the device enters to measuring mode automatically.

CLEANING INSTRUCTIONS:

1. Clean the pacifier thermometer by wiping it with soap and warm water or rubbing alcohol (70% isopropyl alcohol).
2. Wipe dry with soft cloth.
3. Do not boil or clean in a dishwasher. This will cause the thermometer to no longer function and will void the warranty.
4. Do not wash the unit with any thinner or chemical solvent.

SPECIFICATIONS:

FEATURES	DESCRIPTION
Measuring range	90°F ~109.9°F (32°C ~42.9°C)
Resolution	0.1°F (0.1°C)
Accuracy	Comply with ASTM E1112-00 and EN12470-3 requirement
Measurement time	Oral: 1.5 min. (for reference only, it differs from people to people)
Power source	One 3.0V button cell battery (CR2032X1)
Back light	1. Green light is on when the reading is below 99.5°F (37.5°C). 2. Red light is on when the reading is above 99.5°F (37.5°C).
Beeper function	a. 1 beep sound is made when power ON or OFF. b. 10 beeps are made after the measuring completed.
Memory function	Turn on the unit it automatically display the last reading for 2 seconds after display check.
Auto Power-Off function	Around 3 and half min. after the measurement completed.
Out-of-range indication	Displays "Lo" when the measurement is lower than 89.6°F (32°C) and "Hi" when higher than 109.2°F (42.9°C).
Low battery indication	When the battery voltage drops to the low level the symbol "■□" will appear on LCD.
Waterproof function	Yes
Operating temperature	50~104°F (10~40°C), 15~95%R.H. noncondensation
Storage temperature	-4°F to +140°F (-20°C to +60°C)
Self-test function	It runs automatically after power on. Display "Err" when system malfunction.

WARRANTY

We grant you 2 years warranty after the date of purchase. This product has been produced with the greatest care according to international quality standards, established in the European Guideline for Medical Products 93/42/EEC, and it was subjected to strict testing before delivery. Should you nevertheless have reasons of complaint, please send the thermometer together with the warranty certificate, filled out, to the service address given on the back.

Any damage caused by improper handling shall not be covered by the warranty. Batteries and packaging are also excluded from the warranty. All other damage claims excluded. A warranty claim must be submitted within the warranty period. Be sure to include: date of purchase, dealer stamp, and name and address of responsible dealer.

Summer Infant, Inc.

6 Blackstone Valley Place, Suite 206
Lincoln, RI 02865 USA 1-800-263-6237

www.summerinfant.com

© 2005 Summer Infant, Inc.

Please retain information for
future reference.

Summer Infant Europe, LTD.

20 Sparrows Herne
Bushey, Hertfordshire
WD23 1FX UK 44 (0)20 8420 4429
Colors and styles may vary.

MADE IN CHINA

6/05