

Quantifiler™ Kits

**Quantifiler™ Human DNA Quantification Kit and
Quantifiler™ Y Human Male DNA Quantification
Kit**

User's Manual

© Copyright 2006, Applied Biosystems. All rights reserved.

For Research, Forensic, or Paternity Use Only. Not for use in diagnostic procedures.

Information in this document is subject to change without notice. Applied Biosystems assumes no responsibility for any errors that may appear in this document.

APPLIED BIOSYSTEMS DISCLAIMS ALL WARRANTIES WITH RESPECT TO THIS DOCUMENT, EXPRESSED OR IMPLIED, INCLUDING BUT NOT LIMITED TO THOSE OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE. IN NO EVENT SHALL APPLIED BIOSYSTEMS BE LIABLE, WHETHER IN CONTRACT, TORT, WARRANTY, OR UNDER ANY STATUTE OR ON ANY OTHER BASIS FOR SPECIAL, INCIDENTAL, INDIRECT, PUNITIVE, MULTIPLE OR CONSEQUENTIAL DAMAGES IN CONNECTION WITH OR ARISING FROM THIS DOCUMENT, INCLUDING BUT NOT LIMITED TO THE USE THEREOF.

NOTICE TO PURCHASER: LIMITED LICENSE

Use of this product is covered by one or more of the following US patents and corresponding patent claims outside the US: 5,079,352, 6,127,155, 5,677,152 (claims 1–23), 5,773,258 (claims 1–10), 5,210,015, 5,487,972, 5,804,375, 5,538,848, 5,723,591, 5,876,930, 6,030,787, and 6,258,569. The purchase of this product includes a limited, non-transferable immunity from suit under the foregoing patent claims for using only this amount of product solely in forensic and paternity testing, including reporting results of purchaser's activities for a fee or other commercial consideration, and also for the purchaser's own internal research. No right under any other patent claims (such as apparatus or system claims in U.S. Patent No. 6,814,934) is conveyed expressly, by implication, or by estoppel. Further information on purchasing licenses may be obtained from the Director of Licensing, Applied Biosystems, 850 Lincoln Centre Drive, Foster City, California 94404, USA.

TRADEMARKS:

Applied Biosystems, AB (Design), ABI PRISM, AmpF~~STR~~STR, Genotyper, GeneScan, Identifiler, MicroAmp, Profiler Plus, QuantiBlot, Quantifiler, and VIC are registered trademarks, and Applera, FAM, MicroAmp, and ROX are trademarks of Applera Corporation or its subsidiaries in the U.S. and/or certain other countries.

AmpliTag Gold and TaqMan are registered trademarks of Roche Molecular Systems, Inc.

All other trademarks are the sole property of their respective owners.

Part Number 4344790 Rev. D
04/2006

Contents

Preface

How to Use This Guide	vii
Safety	ix
How to Obtain More Information	xiii
How to Obtain Support	xiii

Chapter 1 Overview

Product Overview	1-2
Chemistry Overview	1-3
Instrument Overview	1-7
SDS Software Overview	1-8
Real-Time Data Analysis	1-10
Procedural Overview	1-15
Materials and Equipment	1-16

Chapter 2 Software Setup

Section 2.1 7000 SDS Software Setup	2-3
Overview	2-4
Starting the 7000 SDS	2-5
About Plate Documents	2-7
Setting Up a Plate Document	2-10
Setting Up a Plate Document Template	2-22
Section 2.2 7900HT SDS Software Setup	2-25
Overview	2-26
Starting the 7900HT Real-Time PCR System	2-27
About Plate Documents	2-28
Setting Up a Plate Document	2-31
Setting Up a Plate Document Template	2-40

Chapter 3 PCR Amplification

Preparing the DNA Quantification Standard	3-2
Preparing the Reactions	3-5
Running the Reactions	3-7

Chapter 4 Data Analysis and Results

Section 4.1 7000 SDS Data Analysis	4-3
Analyzing the Plate Document	4-3
Viewing Results	4-4
Section 4.2 7900HT SDS Data Analysis	4-7
Analyzing the Plate Document	4-7
Viewing Results	4-8

Chapter 5 Interpretation of Results

Checking Analysis Settings	5-2
Examining the Standard Curve	5-4
Troubleshooting the Standard Curve	5-6
Using the Internal PCR Control System	5-10
Troubleshooting Amplification Plots	5-12
Assessing Quantity	5-16

Chapter 6 Experiments and Results

Overview	6-2
Section 6.1 ABI PRISM® 7000 Sequence Detection System Validation (SDS Software v1.0)	6-3
Precision	6-4
Reproducibility	6-7
Specificity with a Human DNA Panel	6-10
Specificity with a Non-Human Panel	6-11
Specificity with a Bacterial Pools Panel	6-14

Sensitivity	6-16
Stability	6-17
Mixture Studies	6-21
Degraded DNA Studies	6-23
Comparisons with Other Methods	6-27
Comparison with A ₂₆₀ and Quantiblot Kit	6-27
Comparison with A ₂₆₀ and Dye Intercalation	6-31
Assay Background	6-35
 Section 6.2 Applied Biosystems	
7900HT Real-Time PCR System Validation (SDS Software v2.0) . .	6-37
Precision (7900HT SDS)	6-37
Mixture Studies (7900HT SDS)	6-40
Comparisons with Other Methods (7900HT SDS)	6-41
 Section 6.3 Casework Sample Analysis	6-45
 Section 6.4 Applied Biosystems	
7500 Real-Time PCR System Validation (SDS Software v1.2.3) . . .	6-50
Materials and Methods	6-50
Experimental Setup	6-51
Data Collection	6-54
Data Analysis	6-54
Discussion	6-63
Conclusion	6-64
 Section 6.5 ABI PRISM® 7000 Sequence Detection System Validation	
(SDS Software v1.2.3)	6-66
Materials and Methods	6-66
Experimental Setup	6-67
Data Collection	6-69
Data Analysis	6-69
Discussion	6-79
Conclusion	6-80

Bibliography

Index

Preface

This preface contains:

How to Use This Guide	vii
Safety.	ix
How to Obtain More Information	xiv
How to Obtain Support	xiv

How to Use This Guide

Purpose of This Guide	The <i>Quantifiler® Kits User's Manual</i> provides information about and instructions for using the Quantifiler® Human DNA Quantification Kit and the Quantifiler® Y Human Male DNA Quantification Kit.
------------------------------	--

Text Conventions	<p>This guide uses the following conventions:</p> <ul style="list-style-type: none">• Bold indicates user action. For example: Type 0, then press Enter for each of the remaining fields.• <i>Italic</i> text indicates new or important words and is also used for emphasis. For example: Before analyzing, <i>always</i> prepare fresh matrix.• A right arrow bracket (>) separates successive commands you select from a drop-down or shortcut menu. For example: Select File > Open > Spot Set. Right-click the sample row, then select View Filter > View All Runs.
-------------------------	---

User Attention Words

Two user attention words appear in Applied Biosystems user documentation. Each word implies a particular level of observation or action as described below:

Note: Provides information that may be of interest or help but is not critical to the use of the product.

IMPORTANT! Provides information that is necessary for proper instrument operation, accurate chemistry kit use, or safe use of a chemical.

Examples of the user attention words appear below:

Note: The size of the column affects the run time.

Note: The Calibrate function is also available in the Control Console.

IMPORTANT! To verify your client connection to the database, you need a valid Oracle user ID and password.

IMPORTANT! You must create a separate Sample Entry Spreadsheet for each 96-well plate.

Safety Alert Words

Safety alert words also appear in user documentation. For more information, see [“Safety Alert Words”](#) on [page ix](#).

Safety

Safety Alert Words

Four safety alert words appear in Applied Biosystems user documentation at points in the document where you need to be aware of relevant hazards. Each alert word—**IMPORTANT**, **CAUTION**, **WARNING**, **DANGER**—implies a particular level of observation or action, as defined below:

IMPORTANT! – Indicates information that is necessary for proper instrument operation, accurate chemistry kit use, or safe use of a chemical.



CAUTION – Indicates a potentially hazardous situation that, if not avoided, may result in minor or moderate injury. It may also be used to alert against unsafe practices.



WARNING – Indicates a potentially hazardous situation that, if not avoided, could result in death or serious injury.



DANGER – Indicates an imminently hazardous situation that, if not avoided, will result in death or serious injury. This signal word is to be limited to the most extreme situations.

Chemical Hazard Warning



WARNING CHEMICAL HAZARD. Some of the chemicals used with Applied Biosystems instruments and protocols are potentially hazardous and can cause injury, illness, or death.

Chemical Safety Guidelines

To minimize the hazards of chemicals:

- Read and understand the Material Safety Data Sheets (MSDS) provided by the chemical manufacturer before you store, handle, or work with any chemicals or hazardous materials. (See “[About MSDSs.](#)”)
- Minimize contact with chemicals. Wear appropriate personal protective equipment when handling chemicals (for example, safety glasses, gloves, or protective clothing). For additional safety guidelines, consult the MSDS.
- Minimize the inhalation of chemicals. Do not leave chemical containers open. Use only with adequate ventilation (for example, fume hood). For additional safety guidelines, consult the MSDS.
- Check regularly for chemical leaks or spills. If a leak or spill occurs, follow the manufacturer’s cleanup procedures as recommended on the MSDS.
- Comply with all local, state/provincial, or national laws and regulations related to chemical storage, handling, and disposal.

About MSDSs

Chemical manufacturers supply current Material Safety Data Sheets (MSDSs) with shipments of hazardous chemicals to *new* customers. They also provide MSDSs with the first shipment of a hazardous chemical to a customer after an MSDS has been updated. MSDSs provide the safety information you need to store, handle, transport, and dispose of the chemicals safely.

Each time you receive a new MSDS packaged with a hazardous chemical, be sure to replace the appropriate MSDS in your files.

Obtaining MSDSs

The MSDS for any chemical supplied by Applied Biosystems is available to you free 24 hours a day. To obtain MSDSs:

1. Go to <https://docs.appliedbiosystems.com/msdssearch.html>
2. In the Search field of the MSDS Search page:
 - a. Type in the chemical name, part number, or other information that you expect to appear in the MSDS of interest.
 - b. Select the language of your choice.
 - c. Click **Search**.

3. To view, download, or print the document of interest:
 - a. Right-click the document title.
 - b. Select:
 - **Open** – To view the document
 - **Save Target As** – To download a PDF version of the document to a destination that you choose
 - **Print Target** – To print the document
4. To have a copy of an MSDS sent by fax or e-mail, in the Search Results page:
 - a. Select **Fax** or **Email** below the document title.
 - b. Click **RETRIEVE DOCUMENTS** at the end of the document list.
 - c. Enter the required information.
 - d. Click **View/Deliver Selected Documents Now**.

Note: For the MSDSs of chemicals not distributed by Applied Biosystems, contact the chemical manufacturer.

Chemical Waste Hazard



WARNING CHEMICAL WASTE HAZARD. Some wastes produced by the operation of the instrument or system are potentially hazardous and can cause injury, illness, or death.

Chemical Waste Safety Guidelines

To minimize the hazards of chemical waste:

- Read and understand the Material Safety Data Sheets (MSDSs) provided by the manufacturers of the chemicals in the waste container before you store, handle, or dispose of chemical waste.
- Provide primary and secondary waste containers. (A primary waste container holds the immediate waste. A secondary container contains spills or leaks from the primary container. Both containers must be compatible with the waste material and meet federal, state, and local requirements for container storage.)
- Minimize contact with chemicals. Wear appropriate personal protective equipment when handling chemicals (for example, safety glasses, gloves, or protective clothing). For additional safety guidelines, consult the MSDS.

- Minimize the inhalation of chemicals. Do not leave chemical containers open. Use only with adequate ventilation (for example, fume hood). For additional safety guidelines, consult the MSDS.
- Handle chemical wastes in a fume hood.
- After emptying the waste container, seal it with the cap provided.
- Dispose of the contents of the waste tray and waste bottle in accordance with good laboratory practices and local, state/provincial, or national environmental and health regulations.

Waste Disposal

If potentially hazardous waste is generated when you operate the instrument, you must:

- Characterize (by analysis if necessary) the waste generated by the particular applications, reagents, and substrates used in your laboratory.
- Ensure the health and safety of all personnel in your laboratory.
- Ensure that the instrument waste is stored, transferred, transported, and disposed of according to all local, state/provincial, and/or national regulations.

IMPORTANT! Radioactive or biohazardous materials may require special handling, and disposal limitations may apply.

Biological Hazard Safety



WARNING

BIOHAZARD. Biological samples such as tissues, body fluids, and blood of humans and other animals have the potential to transmit infectious diseases. Follow all applicable local, state/provincial, and/or national regulations. Wear appropriate protective eyewear, clothing, and gloves. Read and follow the guidelines in these publications:

- U.S. Department of Health and Human Services guidelines published in *Biosafety in Microbiological and Biomedical Laboratories* (stock no. 017-040-00547-4; <http://bmbi.od.nih.gov>)

- Occupational Safety and Health Standards, Bloodborne Pathogens (29 CFR§1910.1030; http://www.access.gpo.gov/nara/cfr/waisidx_01/29cfr1910a_01.html).

Additional information about biohazard guidelines is available at:
<http://www.cdc.gov>

How to Obtain More Information

Related Documentation

- *ABI PRISM® 7000 Sequence Detection System User Guide* – Describes the 7000 SDS hardware and software and provides information on preparing, maintaining, and troubleshooting the system.
- *ABI PRISM® 7900HT Sequence Detection System User Guide* – Describes the 7900HT SDS hardware and software and provides information on preparing, maintaining, and troubleshooting the system.

Note: For additional documentation, see “How to Obtain Support” below.

Send Us Your Comments

Applied Biosystems welcomes your comments and suggestions for improving its user documents. You can e-mail your comments to:

techpubs@appliedbiosystems.com

How to Obtain Support

For the latest services and support information for all locations, go to <http://www.appliedbiosystems.com>, then click the link for **Support**.

At the Support page, you can:

- Obtain worldwide telephone and fax numbers to contact Applied Biosystems Technical Support and Sales facilities
- Search through frequently asked questions (FAQs)
- Submit a question directly to Technical Support
- Order Applied Biosystems user documents, MSDSs, certificates of analysis, and other related documents
- Download PDF documents
- Obtain information about customer training
- Download software updates and patches

This chapter covers:

Product Overview	1-2
Chemistry Overview	1-3
Instrument Overview	1-7
SDS Software Overview	1-8
Real-Time Data Analysis	1-10
Procedural Overview	1-15
Materials and Equipment	1-16

Product Overview

Purpose The Quantifiler® Human DNA Quantification Kit (PN 4343895) and the Quantifiler® Y Human Male DNA Quantification Kit (PN 4343906) are designed to quantify the total amount of amplifiable human (and higher primate) DNA or human male DNA in a sample. The results from using the kits can aid in determining:

- If sufficient human DNA or human male DNA is present to proceed with short tandem repeat (STR) analysis
- How much sample to use in STR analysis applications

Product Description The Quantifiler kits contain all the necessary reagents for the amplification, detection, and quantification of a human-specific DNA target or a human male-specific DNA target.

The reagents are designed and optimized for use with the following Applied Biosystems instruments and software:

- ABI PRISM® 7000 Sequence Detection System and SDS Software v1.0
- ABI PRISM® 7900HT Sequence Detection System (no automation module) and SDS Software v2.0.

See [Chapter 6, “Experiments and Results,”](#) for validation studies performed using the Applied Biosystems 7500 Real-Time PCR System with SDS Software v1.2.3 and the ABI PRISM® 7000 Sequence Detection System with SDS Software v1.2.3.

Chemistry Overview

Assay Overview The DNA quantification assay combines two 5' nuclease assays:

- A target-specific (human DNA or human male DNA) assay
- An internal PCR control (IPC) assay

Target-Specific Assay Components The target-specific assay consists of:

- Two primers for amplifying human DNA or human male DNA
- One TaqMan[®] MGB probe labeled with FAM[™] dye for detecting the amplified sequence

About the Targets [Table 1-1](#) provides information about the targets of PCR amplification in the Quantifiler Human kit and the Quantifiler Y kit.

Table 1-1 Targets of Quantifiler kits

Kit	Gene Target	Location	Amplicon Length	Region Amplified	Ploidy
Quantifiler Human kit	Human telomerase reverse transcriptase gene (hTERT)	5p15.33	62 bases	Nontranslated region (intron)	Diploid ^a
Quantifiler Y kit	Sex-determining region Y gene (SRY)	Yp11.3	64 bases	Nontranslated region	Haploid ^a

a. Single-copy target.

IPC Assay Components The IPC assay consists of:

- IPC template DNA (a synthetic sequence not found in nature)
- Two primers for amplifying the IPC template DNA
- One TaqMan[®] MGB probe labeled with VIC[®] dye for detecting the amplified IPC DNA

About the Probes

The TaqMan MGB probes contain:

- A reporter dye (FAM™ dye or VIC® dye) linked to the 5' end of the probe
- A minor groove binder (MGB) at the 3' end of the probe
This modification increases the melting temperature (T_m) without increasing probe length (Afonina *et al.*, 1997; Kuttyavin *et al.*, 1997), which allows the design of shorter probes.
- A nonfluorescent quencher (NFQ) at the 3' end of the probe
Because the quencher does not fluoresce, Applied Biosystems sequence detection systems can measure reporter dye contributions more accurately.

5' Nuclease Assay Process

The 5' nuclease assay process (Figures 1-1 through 1-5) takes place during PCR amplification. This process occurs in every cycle and does not interfere with the exponential accumulation of product.

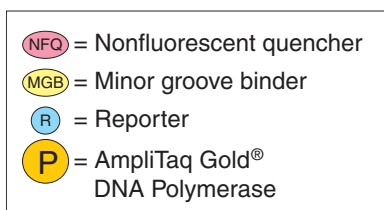


Figure 1-1 Legend for 5' nuclease assay process figures

During PCR, the TaqMan MGB probe anneals specifically to a complementary sequence between the forward and reverse primer sites (Figure 1-2).

When the probe is intact (Figures 1-2 and 1-3), the proximity of the reporter dye to the quencher dye results in suppression of the reporter fluorescence primarily by Förster-type energy transfer (Förster, 1948; Lakowicz, 1983).

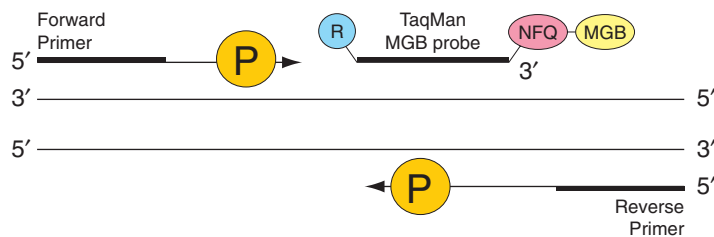


Figure 1-2 Polymerization

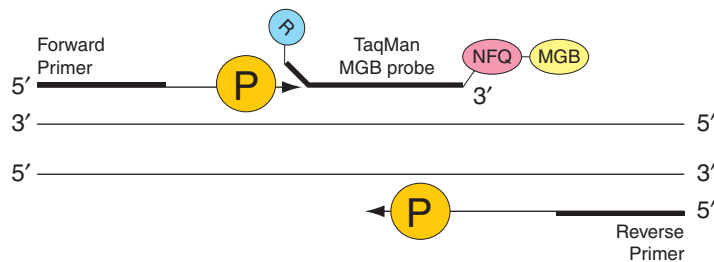


Figure 1-3 Strand displacement

AmpliTaq Gold[®] DNA polymerase cleaves only probes that are hybridized to the target (Figure 1-4). Cleavage separates the reporter dye from the quencher dye, which results in increased fluorescence by the reporter. The increase in fluorescence signal occurs only if the target sequence is complementary to the probe and is amplified during PCR. Because of these requirements, nonspecific amplification is not detected.

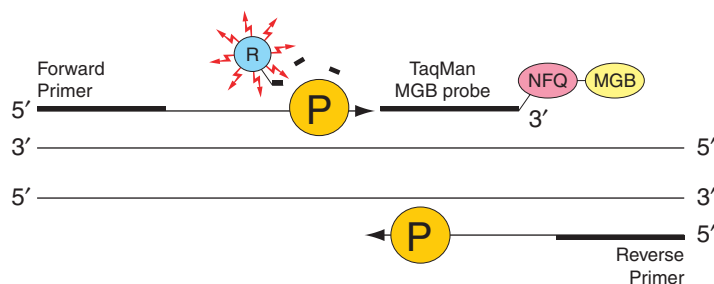


Figure 1-4 Cleavage

Polymerization of the strand continues, but because the 3' end of the probe is blocked, there is no extension of the probe during PCR (Figure 1-5).

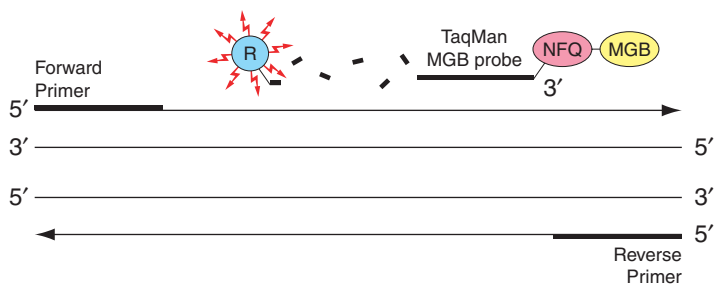


Figure 1-5 Completion of polymerization

Instrument Overview

Fluorescence Detection

Detection on the ABI PRISM 7000 Sequence Detection System

1. A tungsten-halogen lamp directs light to each well on the reaction plate. The light passes through the ABI PRISM™ Optical Adhesive Cover and excites the fluorescent dyes in each well of the plate.
2. A system of lenses, filters, and a dichroic mirror focuses the fluorescence emission into a charge-coupled device (CCD) camera.
3. Based on wavelength, the filters separate the light into a predictably spaced pattern across the CCD camera.
4. During the run, the CCD camera detects the fluorescence emission between 500 nm and 660 nm from each well.
5. The SDS software obtains the fluorescence emission data from the CCD camera and applies data analysis algorithms.

Detection on the ABI PRISM 7900HT Sequence Detection System

1. An argon ion laser directs light to each well on the microplate. The light passes through the ABI PRISM Optical Adhesive Cover and excites the fluorescent dyes in each well of the plate.
2. A system of lenses, filters, and a dichroic mirror focuses the fluorescence emission into a grating.
3. Based on wavelength, the grating separates the light into a predictably spaced pattern across the CCD camera.
4. During the run, the CCD camera detects the fluorescence emission between 500 nm and 660 nm from each well.
5. The SDS software obtains the fluorescence emission data from the CCD camera and applies data analysis algorithms.

SDS Software Overview

This section describes how the SDS software analyzes raw run data from real-time runs. Raw data consists of the spectral data between 500 nm to 660 nm collected by the SDS software during a sequence detection run.

Composite Spectrum

Figure 1-6 shows a composite fluorescence spectrum from a single well containing the passive reference, one probe labeled with FAM[™] dye and a nonfluorescent quencher, and one probe labeled with VIC[®] dye and a nonfluorescent quencher. The example shows how the overlapping component dye spectra contribute to the composite spectrum.

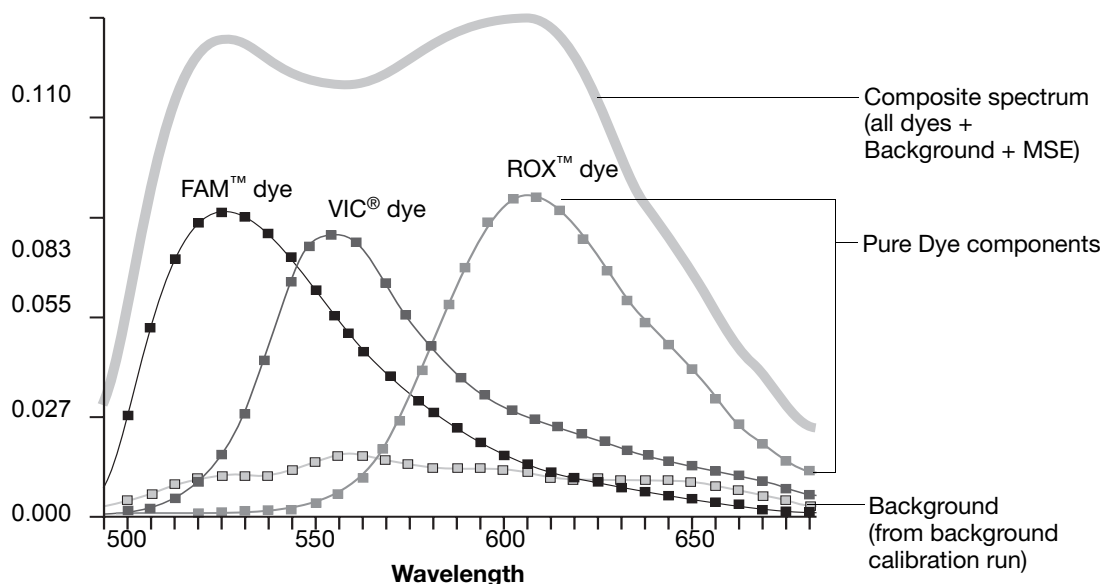


Figure 1-6 Example of a composite spectrum

Processing Multicomponent Data

During the multicomponent transformation, the SDS software uses algorithms to determine the contribution of each dye:

- An algorithm removes the background component stored in the background calibration file to eliminate the contribution of background fluorescence in the raw data.
- The software uses the extracted pure dye standards to express the composite spectrum in terms of the pure dye components.
- Then, an algorithm applies matrix calculations to determine the contributions of each component dye to the composite spectrum.

The software uses the pure dye spectra, generated as part of instrument calibration, to solve for coefficients a, b, and c and to calculate the mean standard error (MSE) in the following equation:

$$\text{Measured spectrum} = a(\text{FAM}) + b(\text{VIC}) + c(\text{ROX}) + d(\text{Background}) + \text{MSE}$$

where coefficients a, b, and c represent the contribution of each dye to the composite spectrum. The MSE value indicates how closely the collective multicomponent spectrum conforms to the raw spectra.

Note: The example equation above assumes that pure dye components exist for FAM dye, VIC dye, and ROX dye and for the instrument background.

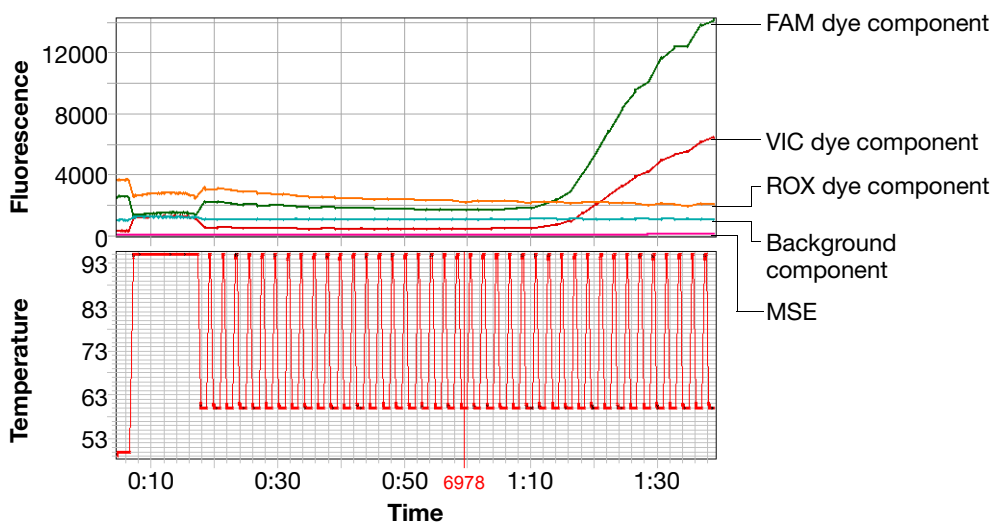


Figure 1-7 Typical component contributions in a multiplex reaction

Normalization of Reporter Signals

The SDS software displays cycle-by-cycle changes in normalized reporter signal (R_n). The SDS software normalizes each reporter signal by dividing it by the fluorescent signal of the passive reference dye. Because the passive reference is one component of the PCR master mix, it is present at the same concentration in all wells of the reaction plate. By normalizing the reporter signal using the passive reference, the software can account for minor variations in signal caused by pipetting inaccuracies and make better well-to-well comparisons of reporter signal.

Real-Time Data Analysis

The 7000 SDS and the 7900HT SDS can be used to determine the relative quantity of a target nucleic acid sequence in a sample by analyzing the cycle-to-cycle change in fluorescent signal as a result of amplification during a PCR ([Figure 1-8](#)).

Amplification Plot Example

When using TaqMan probes with the 7000 SDS or 7900HT SDS, the fluorescent signal (or normalized reporter, R_n) increases as the amount of specific amplified product increases. [Figure 1-8](#) shows the amplification of PCR product in a plot of R_n vs. cycle number during PCR. This amplification plot contains three distinct phases that characterize the progression of the PCR.

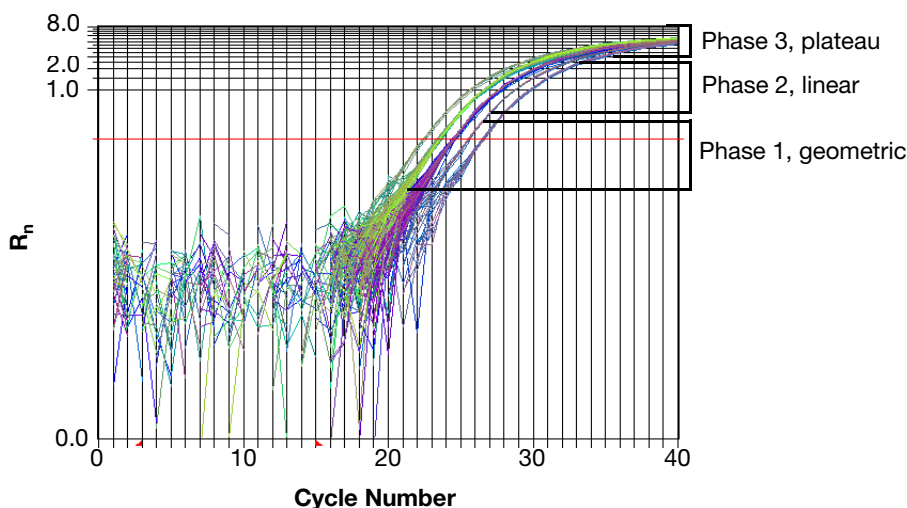


Figure 1-8 Phases of PCR amplification

Phases of Amplification

Initially, R_n appears as a flat line because the fluorescent signal is below the detection limit of the sequence detector.

Phase 1: Geometric (Exponential)

Signal is detected and increases in direct proportion to the increase of PCR product. As PCR product continues to increase, the ratio of AmpliTaq Gold polymerase to PCR product decreases.

During the geometric phase, amplification is characterized by a high and constant efficiency. It occurs between the first detectable rise in fluorescence and before the beginning of the linear phase. During the geometric phase, a plot of DNA concentration versus cycle number on a log scale should approximate a straight line with a slope.

Typically, the SDS system is sufficiently sensitive to detect at least 3 cycles in the geometric phase, assuming reasonably optimized PCR conditions. When the template concentration reaches 10^{-8} M, PCR product stops accumulating exponentially.

Phase 2: Linear

During the linear phase, the slope of the amplification plot decreases steadily. At this point, one or more components of the PCR has decreased below a critical concentration, and the amplification efficiency begins to decrease. This phase is termed linear because amplification approximates an arithmetic progression, rather than a geometric increase. Because the amplification efficiency is continually decreasing during the linear phase, it exhibits low precision.

Phase 3: Plateau

The amplification plot achieves the plateau phase when the PCR stops, the R_n signal remains relatively constant, and the template concentration reaches a plateau at about 10^{-7} M (Martens and Naes, 1989).

Relationship of Amplified PCR Product to Initial Template Concentration

Because of the progressive cleavage of TaqMan fluorescent probes during the PCR, as the concentration of amplified product increases in a sample, so does the R_n value. The following equation describes the relationship of amplified PCR product to initial template during the geometric phase:

$$N_c = N(1 + E)^c$$

where N_c is the concentration of amplified product at any cycle, N is the initial concentration of target template, E is the efficiency of the system, and c is the cycle number.

For example, with the dilutions of RNase P target in the TaqMan[®] RNase P Instrument Verification Plate, the ratio of template concentration to detectable signal is preserved in the geometric phase for all dilutions (Figure 1-9). As the rate of amplification approaches a plateau, the amount of product is no longer proportional to the initial number of template copies.

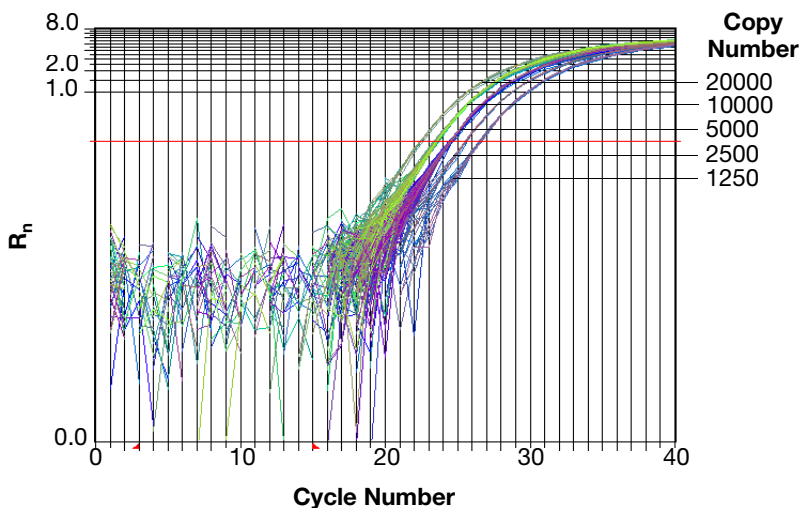


Figure 1-9 Amplification plot from a real-time run of an RNase P Instrument Verification Plate

About the Threshold

The SDS software uses a threshold setting to define the level of detectable fluorescence. Based on the number of cycles required to reach the threshold, the SDS software can compare test samples quantitatively: A sample with a higher starting template copy number reaches the threshold earlier than a sample with a lower starting template copy number.

About the Threshold Cycle

The threshold cycle (C_T) for a specified amplification plot occurs when the fluorescent signal increases beyond the value of the threshold setting. The C_T value depends on:

- Starting template copy number
- Efficiency of DNA amplification by the PCR system

How C_T Values Are Determined

To determine the C_T value, the SDS software uses the R_n values collected from a predefined range of PCR cycles called the baseline (the default baseline occurs between cycles 6 and 15 on the 7000 SDS and between cycles 3 and 15 on the 7900HT SDS):

1. The software generates a baseline-subtracted amplification plot of ΔR_n versus cycle number.
2. An algorithm defines the cycle where the ΔR_n value crosses the threshold setting (the default threshold setting is 0.2) as the threshold cycle (C_T).

**Relationship of
Threshold Cycles
to Initial Template
Amount**

The following equation describes the exponential amplification of the PCR:

$$X_n = X_m(1 + E_X)^{n-m}$$

where:

X_n = number of target molecules at cycle n (so that $n > m$)

X_m = number of target molecules at cycle m

E_X = efficiency of target amplification (between 0 and 1)

$n - m$ = number of cycles elapsed between cycle m and cycle n

Amplicons designed and optimized according to Applied Biosystems guidelines (amplicon size <150 bp) have amplification efficiencies that approach 100%. Therefore $E_X = 1$ so that:

$$\begin{aligned} X_n &= X_m(1 + 1)^{n-m} \\ &= X_m(2)^{n-m} \end{aligned}$$

To define the significance in amplified product of one thermal cycle, set $n - m = 1$ so that:

$$\begin{aligned} X_n &= X_m(2)^1 \\ &= 2X_m \end{aligned}$$

Therefore, each cycle in the PCR reaction corresponds to a two-fold increase in product. Likewise, a difference in C_T values of 1 equates to a two-fold difference in initial template amount.

Procedural Overview

Use of the Quantifiler kits involves the following workflow:

Software Setup



PCR Amplification



Data Analysis



Interpretation of Results

Materials and Equipment

Kit Contents and Storage Each Quantifiler® kit contains materials sufficient to perform 400 reactions at a 25-µL reaction volume.

Table 1-2 Quantifiler kit contents

Reagent	Contents	Quantity	Storage
Quantifiler® Human Primer Mix or Quantifiler® Y Human Male Primer Mix	<ul style="list-style-type: none"> Forward and reverse primers to amplify human DNA or human male DNA target Probe to detect human DNA or human male DNA target IPC system primers, template, and probe 	3 tubes, 1.4 mL each	–15 to –25 °C
Quantifiler Human DNA Standard	200 ng/µL purified DNA standard	1 tube, 120 µL	–15 to –25 °C
Quantifiler PCR Reaction Mix	AmpliTaq Gold® DNA Polymerase, dNTPs with dUTP, Passive Reference, and optimized buffer components	1 tube, 5 mL	2 to 8 °C

Additional Storage Guidelines For Primer Mixes

Follow the additional guidelines for storing the primer mixes:

- Minimize freeze-thaw cycles.
- Keep protected from direct exposure to light. Excessive exposure to light may affect the fluorescent probes.

Equipment and Materials Not Included

Tables 1-3 through 1-5 list required and optional equipment and materials not supplied with the Quantifiler kits. Unless otherwise noted, many of the items are available from major laboratory suppliers (MLS).

Table 1-3 Equipment

Equipment	Source
Applied Biosystems 7900HT Real-Time PCR System (no automation)	Contact your local Applied Biosystems sales representative.
ABI PRISM® 7000 Sequence Detection System	
Tabletop centrifuge with 96-well plate adapters (optional)	MLS

Table 1-4 User-supplied materials

Material	Source
Quantifiler® Human DNA Quantification Kit	Applied Biosystems (PN 4343895)
Quantifiler® Y Human Male DNA Quantification Kit	Applied Biosystems (PN 4343906)
Glycogen, 20 mg (1 mL)	Roche Applied Science (PN 901 393)
High-Throughput Setup	
96-Well Optical Reaction Plates	Applied Biosystems (PN 4306737)
Optical Adhesive Covers Starter Kit (20 covers, 1 compression pad, 1 applicator)	Applied Biosystems (PN 4313663)
Optical Adhesive Covers (100 covers)	Applied Biosystems (PN 4311971)
MicroAmp™ Splash Free Support Base	Applied Biosystems (PN 4312063)

Table 1-4 User-supplied materials (*continued*)

Material	Source
Mid-to-Low-Throughput Setup	
MicroAmp® Optical Tubes (8 tubes/strip, 125 strips)	Applied Biosystems (PN 4316567)
MicroAmp™ 96-Well Tray/Retainer Set	Applied Biosystems (PN 403081)
Optical Caps (8 caps/strip, 300 strips)	Applied Biosystems (PN 4323032)
Compression pad from Optical Adhesive Covers Starter Kit Note: Not necessary if using Optical Caps	Applied Biosystems (PN 4313663)

Table 1-5 Documents

Document	Applied Biosystems Part Number
<i>ABI PRISM® 7000 Sequence Detection System User Guide</i>	4330228
<i>ABI PRISM® 7900HT Sequence Detection System User Guide</i>	4317596

Chapter 2

Software Setup

This chapter covers:

Section 2.1 7000 SDS Software Setup	2-3
Overview	2-4
Starting the 7000 SDS	2-5
About Plate Documents	2-7
Setting Up a Plate Document	2-10
Setting Up a Plate Document Template	2-22
Section 2.2 7900HT SDS Software Setup	2-25
Overview	2-26
Starting the 7900HT Real-Time PCR System	2-27
About Plate Documents	2-28
Setting Up a Plate Document	2-31
Setting Up a Plate Document Template	2-40

Section 2.1 7000 SDS Software Setup

This section covers:

Overview	2-4
Starting the 7000 SDS.....	2-5
About Plate Documents.....	2-7
Setting Up a Plate Document	2-10
Setting Up a Plate Document Template	2-22

Overview

Purpose During software setup, you start up the ABI PRISM® 7000 Sequence Detection System (7000 SDS) and set up a plate document for DNA quantification using the Quantifiler® kits.


Configuration The Quantifiler® kits are supported using the 7000 SDS and Sequence Detection Systems (SDS) Software v1.0 for real-time data collection and analysis.

Starting the 7000 SDS

Overview Starting the 7000 SDS involves:

1. [Starting the Computer](#)
2. [Powering On the Instrument \(page 2-6\)](#)
3. [Starting SDS Software \(page 2-6\)](#)

Starting the Computer

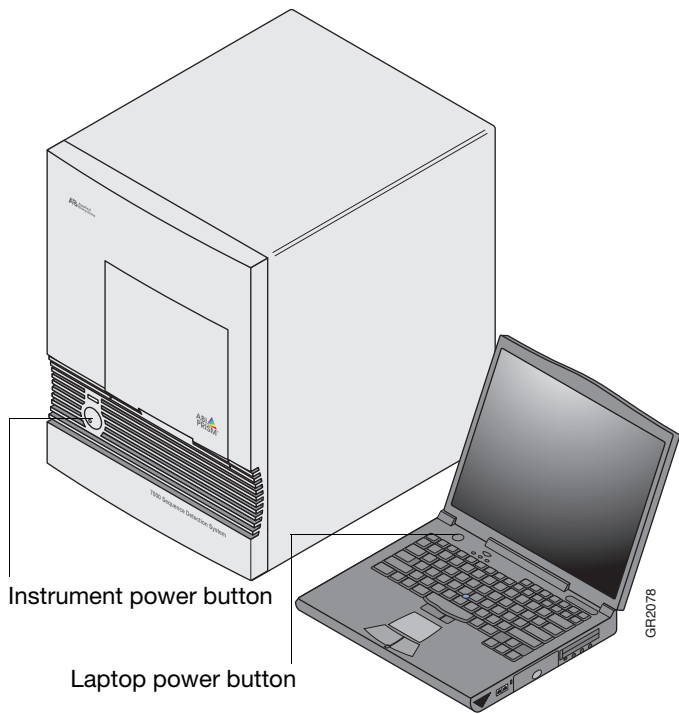1.	If you are using the laptop computer, open it by pushing in the front, center button, holding it, and lifting up the lid.
2.	Press the power button on the computer.  <p>The diagram illustrates the physical components involved in starting the computer. On the left is the 7000 SDS instrument, a large, light-colored rectangular unit with a front panel featuring a power button and a label that reads '7000 SDS Instrument'. To the right is a laptop computer, shown open. A line points from the text 'Instrument power button' to the power button on the instrument's front panel. Another line points from the text 'Laptop power button' to the power button on the laptop's keyboard area. The laptop is labeled 'GR2078' on its side.</p>
3.	In the Enter User name field of the login window, type your name or the user name associated with the computer.
4.	If required, type your password in the Password field.

Powering On the Instrument

Note: Wait for the computer to finish starting up before powering on the 7000 instrument.



WARNING

PHYSICAL INJURY HAZARD. Moving parts can crush and cut. Keep hands clear of moving parts while operating. Disconnect power before servicing the instrument.


Press the power button on the lower left front of the instrument.



Power button

Starting SDS Software

Select **Start > ABI Prism 7000 > ABI Prism 7000 SDS Software**.

The software attempts to initialize the instrument and displays a message in the status bar for a few seconds. Then the computer attempts to establish communication with the 7000 instrument. If the connection is successful, the software displays  in the status bar.

About Plate Documents

How Plate Documents Are Used

Running a reaction plate on the 7000 SDS requires creating and setting up a plate document using the SDS software. A plate document is a representation of the arrangement of samples (standards and unknowns) and detectors on the reaction plate. The SDS software uses the plate document to:

- Coordinate the instrument operation, such as thermal cycling and data collection
- Organize and store the data gathered during the run
- Analyze the data from the run

Plate Document Types

You can use the SDS software to create two types of plate document files.

Plate Document Type	File Extension	Description
SDS document	*.sds	Primary file to use when performing a run. Required for all experiments.
SDS template	*.sdt	File that already contains run parameters that are commonly used in plate documents, such as detectors, thermal cycler conditions, and so on. Streamlines the creation of the SDS document (*.sds) file.

Example Plate Document Setup

You can arrange the reactions in any well of the reaction plate, but you need to set up the plate document so that it corresponds exactly to the arrangement of the standards and unknown samples in the wells of the reaction plate.

Figure 2-1 shows one example of arranging reactions when running two Quantifiler kits on one 96-well reaction plate:

- Wells A1 through D12 (gray) correspond to reactions using the Quantifiler Human kit
- Wells E1 through H12 (white) correspond to reactions using the Quantifiler Y kit

Note: For each Quantifiler kit assay, there are eight DNA quantification standards and two reactions for each standard. See [“Preparing the DNA Quantification Standard” on page 3-2](#) for more information about the DNA quantification standards.

	1	2	3	4	5	6	7	8	9	10	11	12
A	Std 1	Std 1	Std 2	Std 2	Std 3	Std 3	Std 4	Std 4	Std 5	Std 5	Std 6	Std 6
B	Std 7	Std 7	Std 8	Std 8	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN
C	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN
D	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN
E	Std 1	Std 1	Std 2	Std 2	Std 3	Std 3	Std 4	Std 4	Std 5	Std 5	Std 6	Std 6
F	Std 7	Std 7	Std 8	Std 8	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN
G	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN
H	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN

Figure 2-1 Example plate setup of reactions with two kits

Figure 2-2 shows another example of arranging reactions when running two Quantifiler kits on one 96-well reaction plate if you are using repeat pipettors:

- Wells A1 through D6 (gray) correspond to reactions using the Quantifiler Human kit
- Wells A7 through H12 (white) correspond to reactions using the Quantifiler Y kit

Note: For each Quantifiler kit assay, there are eight DNA quantification standards and two reactions for each standard. See [“Preparing the DNA Quantification Standard” on page 3-2](#) for more information about the DNA quantification standards.

	1	2	3	4	5	6	7	8	9	10	11	12
A	Std 1	Std 1	UNKN	UNKN	UNKN	UNKN	Std 1	Std 1	UNKN	UNKN	UNKN	UNKN
B	Std 2	Std 2	UNKN	UNKN	UNKN	UNKN	Std 2	Std 2	UNKN	UNKN	UNKN	UNKN
C	Std 3	Std 3	UNKN	UNKN	UNKN	UNKN	Std 3	Std 3	UNKN	UNKN	UNKN	UNKN
D	Std 4	Std 4	UNKN	UNKN	UNKN	UNKN	Std 4	Std 4	UNKN	UNKN	UNKN	UNKN
E	Std 5	Std 5	UNKN	UNKN	UNKN	UNKN	Std 5	Std 5	UNKN	UNKN	UNKN	UNKN
F	Std 6	Std 6	UNKN	UNKN	UNKN	UNKN	Std 6	Std 6	UNKN	UNKN	UNKN	UNKN
G	Std 7	Std 7	UNKN	UNKN	UNKN	UNKN	Std 7	Std 7	UNKN	UNKN	UNKN	UNKN
H	Std 8	Std 8	UNKN	UNKN	UNKN	UNKN	Std 8	Std 8	UNKN	UNKN	UNKN	UNKN

Figure 2-2 Example plate setup of reactions using repeat pipettors

Setting Up a Plate Document


7000 SDS Software Setup

Overview Setting up a plate document to run Quantifiler kit assays involves:

- 1. [Creating a Blank Plate Document](#) (page 2-10)
- 2. [Creating Detectors](#) (the first time only, page 2-11)
- 3. [Adding Detectors to the Plate Document](#) (page 2-14)
- 4. [Applying Detectors for Standards](#) (page 2-15)
- 5. [Applying Detectors for Unknown Samples](#) (page 2-17)
- 6. [Adding Sample Names for Unknown Samples](#) (page 2-18)
- 7. [Setting Thermal Cycler Conditions](#) (page 2-19)
- 8. [Saving the Plate Document](#) (page 2-21)


Creating a Blank Plate Document

To create a blank plate document:

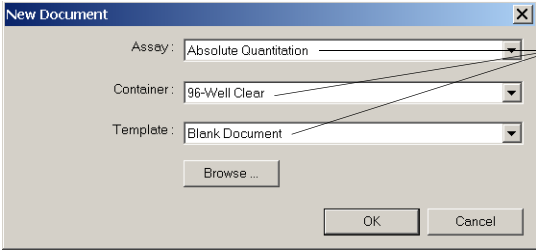1.	If the SDS software is not already started, select Start > Programs > ABI Prism 7000 > ABI Prism 7000 SDS Software .
2.	In the SDS software, select File > New to open the New Document dialog box. <div><p>Default settings</p></div>

To create a blank plate document: *(continued)*

3. Click **OK** to use the default settings and to view a blank plate document:

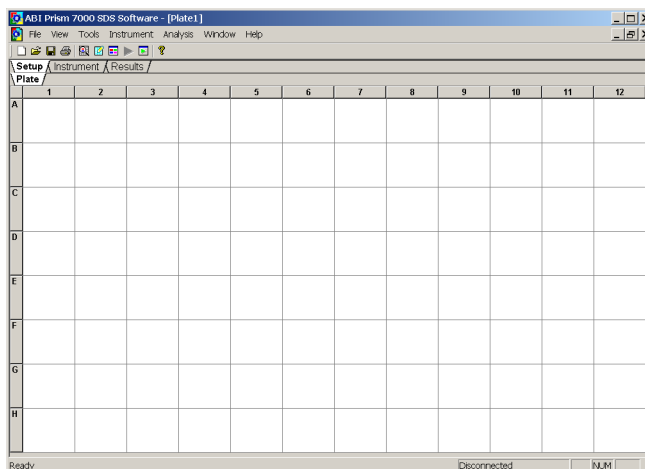

Creating Detectors

Before you set up the plate document, you need to create detectors in the SDS software for running Quantifiler kit assays. After the detectors are created, you do not need to create detectors for subsequent runs of Quantifiler kit assays and you can skip to [“Adding Detectors to the Plate Document” on page 2-14](#).

To create detectors:

1. Select **Tools > Detector Manager**.
2. In the lower left part of the Detector Manager dialog box, select **File > New** to open the New Detector dialog box.

To create detectors: (continued)

3. Create a detector for the Quantifiler Human kit:

The 'New Detector' dialog box is shown with the following fields and values:

- Name: Quantifiler Human
- Description: (empty)
- Reporter Dye: FAM
- Quencher Dye: (none)
- Color: Blue square
- Notes: (empty)

Buttons at the bottom: Create Another, OK, Cancel.

Enter **Quantifiler Human**

Select **FAM**

Make sure **(none)** is selected

Click to select a color

4. Click **Create Another** to add the Quantifiler Human detector and to reset the New Detector dialog box.

5. Create a detector for the Quantifiler Y kit:

The 'New Detector' dialog box is shown with the following fields and values:

- Name: Quantifiler Y
- Description: (empty)
- Reporter Dye: FAM
- Quencher Dye: (none)
- Color: Green square
- Notes: (empty)

Buttons at the bottom: Create Another, OK, Cancel.

Enter **Quantifiler Y**

Select **FAM**

Make sure **(none)** is selected

Click to select a color

6. Click **Create Another** to add the Quantifiler Y detector and to reset the New Detector dialog box.

To create detectors: (continued)

7. Create a detector for the IPC assay:

New Detector

Name: Enter **IPC**

Description:

Reporter Dye: Select **VIC**

Quencher Dye: Make sure **(none)** is selected

Color: Click to select a color

Notes:

8. Click **OK** to add the IPC detector and to return to the Detector Manager dialog box.

Detector Manager

Detector List


Find:

Detector Name	Description	Reporter	Quencher	Color	Notes	Last
Quantifier Human		FAM	(none)			2003/8/20 13:
Quantifier Y		FAM	(none)			2003/8/20 13:
IPC		VIC	(none)			2003/8/20 13:


Adding Detectors to the Plate Document

To add detectors to the plate document:

1. In the SDS software, select **Tools > Detector Manager**. If the detectors for the Quantifiler kits have been created, they are listed in the Detector Manager:

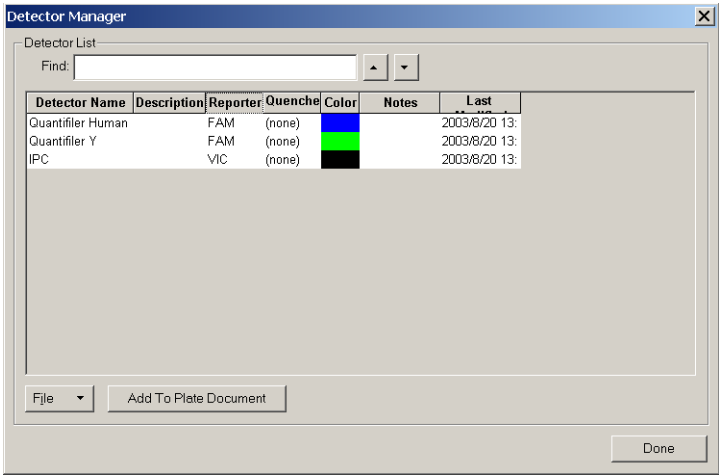

2. In the Detector Manager, select the **Quantifiler Human**, **Quantifiler Y**, and the **IPC** detectors by clicking them while pressing the Ctrl key.

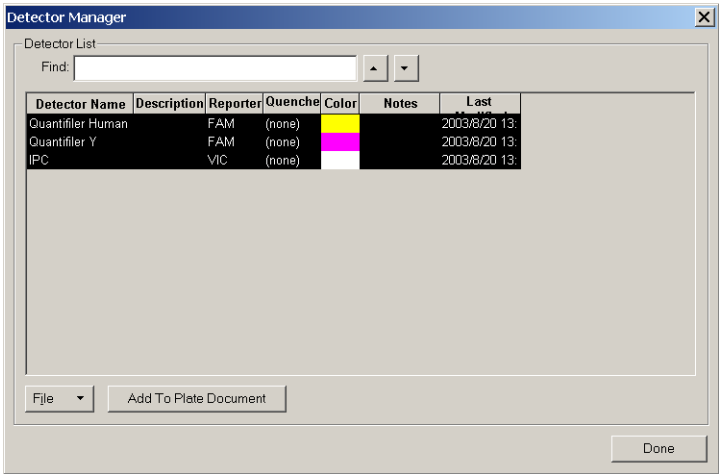

3. Click **Add To Plate Document**.
4. Click **Done** to close the Detector Manager.


Applying Detectors for Standards

You need to apply detectors to the plate document for the wells on the reaction plate that contain DNA quantification standards. Repeat the procedure until you complete applying detector tasks, quantities, and sample names for all quantification standards.

IMPORTANT! Set up detectors for each quantity and for each kit separately. For example, set up detectors for quantification standard 1 for the Quantifiler Human kit first, and then for quantification standard 2 for the Quantifiler Human kit, and so on, until you finish setting up the detectors for all wells containing quantification standards

To apply detectors for quantification standards:


1. Select **View > Well Inspector** to open the dialog box:

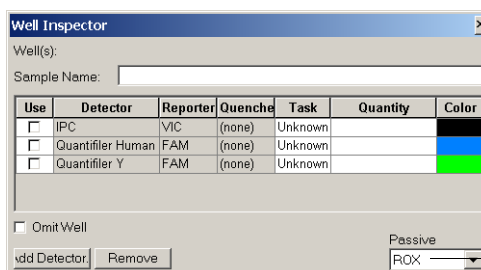

Make sure that **ROX** is selected

Note: The Well Inspector displays the detectors that were added to the plate document.

2. On the Plate tab, select wells that correspond to a specific quantification standard for one kit.

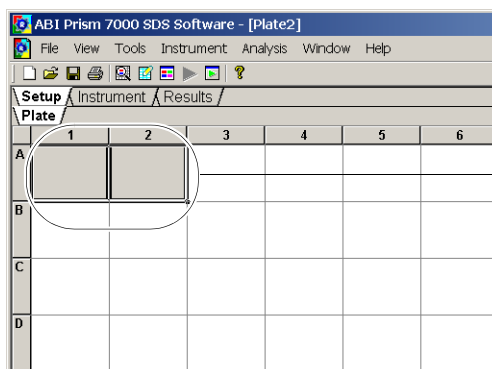

Wells selected

To apply detectors for quantification standards: *(continued)*

3. With the wells selected, go to the Well Inspector and:
 - a. Select the Use boxes for the applicable detectors:
 - IPC
 - Quantifiler Human or Quantifiler Y
 - b. For the Quantifiler Human or Quantifiler Y detector, click **Unknown** in the Task column, then select **Standard** from the drop-down list.
 - c. For the Quantifiler Human or Quantifiler Y detector, select the Quantity field for the appropriate detector and enter the quantity of DNA in the well.

IMPORTANT! Although you do not enter units for Quantity, you must use a consistent unit (for example, ng/ μ L) for all standard quantities. The units used for standard quantities defines the quantification units for analysis results.

Note: Leave the IPC detector Task for standard reactions set to Unknown. Quantity values are not needed for IPC detectors.

- d. Enter the Sample Name (for example, Std. 1, Std. 2, and so on).

For example:

The screenshot shows the 'Well Inspector' dialog box. At the top, it says 'Well(s): A1-A2' and 'Sample Name: Std. 1'. Below this is a table with columns: Use, Detector, Reporter, Quencher, Task, Quantity, and Color. The table has three rows: IPC (checked), Quantifiler Human (checked), and Quantifiler Y (unchecked). The Task for IPC is 'Unknown', and for Quantifiler Human and Y, it is 'Standard'. The Quantity for Quantifiler Human is '50'. The Color column shows black for IPC, blue for Quantifiler Human, and green for Quantifiler Y. Below the table, there is a checkbox for 'Omit Well' and a 'Passive' dropdown menu set to 'ROX'. There are also 'Add Detector' and 'Remove' buttons.

Use	Detector	Reporter	Quencher	Task	Quantity	Color
<input checked="" type="checkbox"/>	IPC	VIC	(none)	Unknown		
<input checked="" type="checkbox"/>	Quantifiler Human	FAM	(none)	Standard	50	
<input type="checkbox"/>	Quantifiler Y	FAM	(none)	Unknown		

Task for IPC set to **Unknown** (default)

Applying Detectors for Unknown Samples

You need to apply detectors to the plate document for the wells on the reaction plate that contain unknown samples.

IMPORTANT! If you run reactions for the Quantifiler Human kit and the Quantifiler Y kit on the same plate, apply detectors for unknown samples for each kit separately.

To apply detectors for unknown samples:

1. On the Plate tab, select the wells that correspond to all unknown samples for one Quantifiler kit.
2. With the well(s) selected, select **View > Well Inspector** and check the Use boxes for the applicable detectors:
 - Quantifiler Human or Quantifiler Y
 - IPC

For example:

Use	Detector	Reporter	Quenche	Task	Quantity	Color
<input checked="" type="checkbox"/>	Quantifiler Human	FAM	(none)	Unknown		Blue
<input type="checkbox"/>	Quantifiler Y	FAM	(none)	Unknown		Green
<input checked="" type="checkbox"/>	IPC	VIC	(none)	Unknown		Black

☐ Omit Well

Add Detector... Remove

Passive
ROX

Make sure that **ROX** is selected

3. If you are running both kits on the reaction plate, repeat [steps 1](#) and [2](#) for the unknown samples for the other kit.
4. Select **View > Well Inspector** to close the Well Inspector.

Adding Sample Names for Unknown Samples

Repeat the procedure to add sample names for all unknown samples.

To add sample names for unknown samples:

1.
- On the Plate tab, select one well containing an unknown sample.
2.
- With the well selected, select **View > Well Inspector** and enter the Sample Name.

For example:

Well Inspector

Well(s): B5

Sample Name: Unknown 1

Use	Detector	Reporter	Quencher	Task	Quantity	Color
<input checked="" type="checkbox"/>	Quantifiler Human	FAM	(none)	Unknown		Blue
<input type="checkbox"/>	Quantifiler Y	FAM	(none)	Unknown		Green
<input checked="" type="checkbox"/>	IPC	VIC	(none)	Unknown		Black

☐ Omit Well

Add Detector

Remove

Passive

ROX

Note: Samples with identical sample names are treated as replicates by the SDS software. Results for replicate reactions are grouped together automatically for data analysis.

Setting Thermal Cycler Conditions

Before running a Quantifiler kit assay, you need to make two changes to the default thermal cycler conditions:

- Thermal profile
- Sample volume

To set thermal cycler conditions:

1. In the plate document, select the **Instrument** tab.

2. Press the **Shift** key and click within the Stage 1 hold step (50 °C for 2 minutes) to select it.

Thermal Profile

Auto Increment

Stage 1

Stage 2

Stage 3

Reps: 1

Reps: 1

Reps: 40

50.0

55.0

55.0

10.00

0.15

60.0

1.00

50.0

2.00

Add Cycle

Add Hold

Add Step

Sample Volume: 50 (uL)

☐ Dissociation Protocol

☒ 9600 Emulation

60.0

C

Press the **Shift** key while you click within the Stage 1 hold step

3. After the hold step is selected, press the **Delete** key.

Thermal Profile

Auto Increment

Stage 1

Stage 2

Stage 3

Reps: 1

Reps: 1

Reps: 40

50.0

55.0

55.0

10.00

0.15

60.0

1.00

50.0

2.00

Add Cycle

Add Hold

Add Step

Sample Volume: 50 (uL)

☐ Dissociation Protocol

☒ 9600 Emulation

60.0

C

Stage 1 hold step is selected


7000 SDS Software Setup

Quantifiler Kits User's Manual

2-19

To set thermal cycler conditions: *(continued)*

4. Make sure that the thermal profile appears as follows:

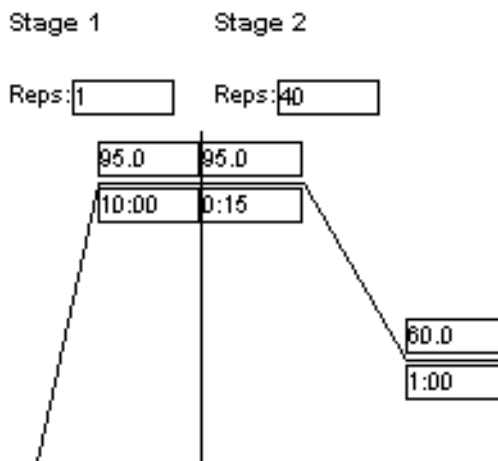

5. Change the Sample Volume to **25 (μL)** and make sure that the 9600 Emulation box is selected.

Note: Selecting the 9600 Emulation box reduces the ramp rate.

Thermal Profile Auto Increment

Stage 1 Stage 2

Reps: 1 Reps: 40

95.0 95.0

10:00 0:15

60.0

1:00

Add Cycle Add Hold Add Step

Sample Volume: 25 (uL)

Dissociation Protocol: 60.0 C

☐ Dissociation Protocol

☒ 9600 Emulation

Set the volume to 25 μL

Make sure that this box is selected

Saving the Plate Document

Before running the reaction plate, save the plate document as an SDS Document (*.sds) file.

Note: To save the plate document as a template, see [“Setting Up a Plate Document Template”](#) on page 2-22.

To save the plate document:

1.	Select File > Save .
2.	Select the location for the plate document.
3.	Enter a file name.
4.	For Save as type, select SDS Documents (*.sds) .
5.	Click Save .

Setting Up a Plate Document Template

Purpose A plate document template reduces the time required to set up a plate document. This section describes how to create an SDS Template Document (*.sdt) for running Quantifiler kit assays.


Template Settings In addition to plate document settings (assay and container), templates can contain:

- Assay-specific detectors
- Well assignments for quantification standards, with detectors, tasks, and quantity
- Well assignments for unknown samples, with detectors and tasks
- Instrument settings: thermal cycler conditions and reaction volume settings

Creating a Plate Document Template

This procedure assumes that you have created the detectors for running reactions using the Quantifiler kits ([page 2-11](#)).

To create a plate document template:

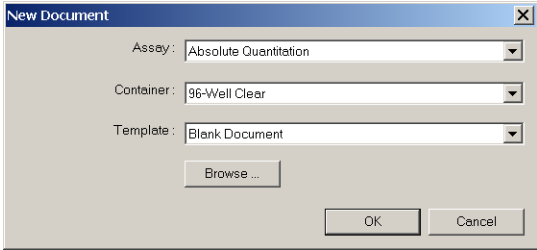1.	If the SDS software is not already started, select Start > Programs > ABI Prism 7000 > ABI Prism 7000 SDS Software .
2.	<p>Select File > New, complete the New Document dialog box, then click OK.</p> 
3.	<p>Apply the desired template settings to the plate document:</p> <ul style="list-style-type: none"> • Add detectors to the plate document (page 2-14) • Apply detectors for standards and for unknown samples (page 2-15 and page 2-17) • Set thermal cycler conditions (page 2-19)


To create a plate document template: (continued)

4. Select **File > Save As** and complete the Save As dialog box:
 - a. For Save as type, select **SDS Templates (*.sdt)**.
 - b. Locate and select the Templates folder within the software folder:

X:Program Files > ABI Prism 7000 > Templates, where *X* is the hard drive on which the SDS software is installed.

Note: Saving the template file in the Templates folder makes the template available in the Template drop-down list of the New Document dialog box (see [step 2](#) in “Creating a Plate Document from a Template” on page 2-24).

- c. For File name, enter a name for the template. For example, enter **Quantifiler Template**:

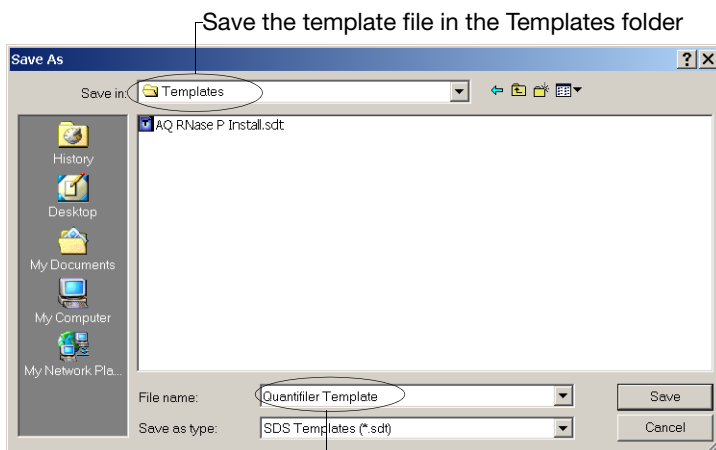

- d. Click **Save**.

Creating a Plate Document from a Template

After you create a template, you can use it to create a plate document.

To create a plate document from a template:

1.	If the SDS software is not already started, select Start > Programs > ABI Prism 7000 > ABI Prism 7000 SDS Software .
2.	<p>Select File > New and in the New Document dialog box and make the following selections:</p> <ul style="list-style-type: none"> • For Assay, select Absolute Quantitation. • For Container, select 96-Well Clear. • For Template, select an appropriate template from the list. <p>Note: If the template is not available in the list, click Browse to locate and select an appropriate template.</p>
3.	<p>Complete the plate document setup:</p> <ul style="list-style-type: none"> • Add detectors to the plate document (page 2-14) • Apply detectors for standards and for unknown samples (page 2-15 and page 2-17) • Set thermal cycler conditions (page 2-19) <p>Note: The tasks that you perform vary according to which settings were defined in the template.</p>
4.	<p>Save the plate document (page 2-21).</p> <p>Note: For Save as type, select SDS Documents (*.sds).</p>

Section 2.2 7900HT SDS Software Setup

This section covers:

Overview	2-26
Starting the 7900HT Real-Time PCR System.	2-27
About Plate Documents.	2-28
Setting Up a Plate Document	2-31
Setting Up a Plate Document Template	2-40

Overview

Purpose During software setup, you start up the Applied Biosystems 7900HT Real-Time PCR System and set up a plate document for DNA quantification using the Quantifiler® kits.

Configuration The Quantifiler® kits are supported using the following configuration of the 7900HT Real-Time PCR System for real-time data collection and analysis:

- 96-well reaction plates
- Manual setup
- Sequence Detection Systems (SDS) software v2.0

Note: Use of the robotic microplate handler and/or 384-well reaction plates is not supported.



Starting the 7900HT Real-Time PCR System

Overview Starting the Applied Biosystems 7900HT Real-Time PCR System involves:

1. Powering on the computer.
2. Powering on the instrument.
3. Starting the SDS software.

Starting the 7900HT System

To start the 7900HT System:

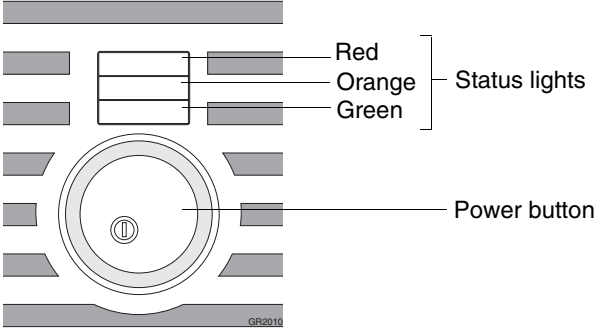1.	Press the power buttons on the computer and on the monitor.
2.	In the login screen, enter the User Name and Password.
3.	<p>Press the power button below the status lights on the front of the instrument.</p>  <p>At startup, the instrument:</p> <ul style="list-style-type: none"> • Emits a high-pitched tone, indicating that the system is initialized • Cycles the status lights (red > orange > green), indicating that the instrument is active
4.	<p>Select Start > Programs > Applied Biosystems > SDS 2.0.</p> <p>At startup, the software attempts to establish communication with the 7900HT instrument. If the connection is successful, the software displays  in the status bar.</p>

About Plate Documents

How Plate Documents Are Used

Running a reaction plate on the 7900HT Real-Time PCR System requires creating and setting up a plate document using the SDS software. A plate document is a representation of the arrangement of samples (standards and unknowns) and reagents on the reaction plate. The SDS software uses the plate document to:

- Coordinate the instrument operation, such as thermal cycling and data collection
- Organize and store the data gathered during the run
- Analyze the data from the run

Plate Document Types

You can use SDS software to create two types of plate document files.

Plate Document Type	File Extension	Description
Single plate document	*.sds	Primary file to use when performing a run. Required for all experiments.
Template plate document	*.sdt	File that already contains run parameters that are commonly used in plate documents, such as detectors, thermal cycler conditions, and so on. Streamlines the creation of the SDS document (*.sds) file.

Example Plate Document Setup

You can arrange the reactions in any well of the reaction plate, but you need to set up the plate document so that it corresponds exactly to the arrangement of the standards and unknown samples in the wells of the reaction plate.

Figure 2-3 shows one example of arranging reactions when running two Quantifiler kit assays on one 96-well plate:

- Wells A1 through D12 (gray) correspond to reactions using the Quantifiler Human kit
- Wells E1 through H12 (white) correspond to reactions using the Quantifiler Y kit

Note: For each Quantifiler kit assay, there are eight DNA quantification standards and two reactions for each standard. See [“Preparing the DNA Quantification Standard” on page 3-2](#) for more information about the DNA quantification standards.

	1	2	3	4	5	6	7	8	9	10	11	12
A	Std 1	Std 1	Std 2	Std 2	Std 3	Std 3	Std 4	Std 4	Std 5	Std 5	Std 6	Std 6
B	Std 7	Std 7	Std 8	Std 8	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN
C	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN
D	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN
E	Std 1	Std 1	Std 2	Std 2	Std 3	Std 3	Std 4	Std 4	Std 5	Std 5	Std 6	Std 6
F	Std 7	Std 7	Std 8	Std 8	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN
G	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN
H	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN	UNKN

Figure 2-3 Example arrangement of reactions with two kits

Figure 2-4 shows another example of arranging reactions when running two Quantifiler kits on one 96-well reaction plate if you are using repeat pipettors:

- Wells A1 through D6 (gray) correspond to reactions using the Quantifiler Human kit
- Wells A7 through H12 (white) correspond to reactions using the Quantifiler Y kit

Note: For each Quantifiler kit assay, there are eight DNA quantification standards and two reactions for each standard. See “Preparing the DNA Quantification Standard” on page 3-2 for more information about the DNA quantification standards.

	1	2	3	4	5	6	7	8	9	10	11	12
A	Std 1	Std 1	UNKN	UNKN	UNKN	UNKN	Std 1	Std 1	UNKN	UNKN	UNKN	UNKN
B	Std 2	Std 2	UNKN	UNKN	UNKN	UNKN	Std 2	Std 2	UNKN	UNKN	UNKN	UNKN
C	Std 3	Std 3	UNKN	UNKN	UNKN	UNKN	Std 3	Std 3	UNKN	UNKN	UNKN	UNKN
D	Std 4	Std 4	UNKN	UNKN	UNKN	UNKN	Std 4	Std 4	UNKN	UNKN	UNKN	UNKN
E	Std 5	Std 5	UNKN	UNKN	UNKN	UNKN	Std 5	Std 5	UNKN	UNKN	UNKN	UNKN
F	Std 6	Std 6	UNKN	UNKN	UNKN	UNKN	Std 6	Std 6	UNKN	UNKN	UNKN	UNKN
G	Std 7	Std 7	UNKN	UNKN	UNKN	UNKN	Std 7	Std 7	UNKN	UNKN	UNKN	UNKN
H	Std 8	Std 8	UNKN	UNKN	UNKN	UNKN	Std 8	Std 8	UNKN	UNKN	UNKN	UNKN

Figure 2-4 Example arrangement of reactions using repeat pipettors

Setting Up a Plate Document

Overview Setting up a plate document involves:

1. [Creating a Blank Plate Document \(page 2-31\)](#)
2. [Creating Detectors \(page 2-32\)](#)
3. [Copying Detectors to the Plate Document \(page 2-34\)](#)
4. [Applying Detectors for Standards \(page 2-35\)](#)
5. [Applying Detectors for Unknown Samples \(page 2-36\)](#)
6. [Adding Sample Names for Unknown Samples \(page 2-37\)](#)
7. [Setting Thermal Cycler Conditions \(page 2-38\)](#)
8. [Saving the Plate Document \(page 2-39\)](#)

Creating a Blank Plate Document


To create a blank plate document:

1. If the SDS software is not already started, select **Start > Programs > Applied Biosystems > SDS 2.0**.
2. Select **File > New**, complete the New Document dialog box, then click **OK**.

**Creating
Detectors**


Before you set up the plate document, you need to create detectors in the SDS software for running Quantifiler kit assays. After the detectors are created, you do not need to create detectors for subsequent runs of Quantifiler kit assays and you can skip to [“Copying Detectors to the Plate Document” on page 2-34](#).

To create detectors:

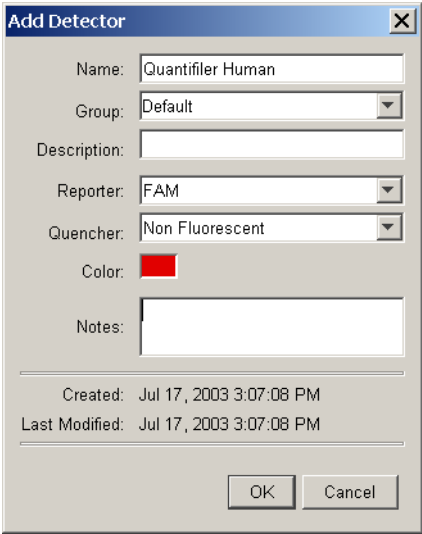- | | |
|----|---|
| 1. | With a new plate document open, select Tools > Detector Manager . |
| 2. | <p>Create a detector for the Quantifiler Human kit:</p> <p>a. In the lower left part of the Detector Manager, click New, then complete the dialog box:</p> <div data-bbox="482 635 904 1171"></div> <p>b. Click OK to return to the Detector Manager.</p> |

To create detectors: (continued)

3. Create a detector for the Quantifiler Y Human Male kit:
 - a. In the Detector Manager, click **New** and complete the dialog box:

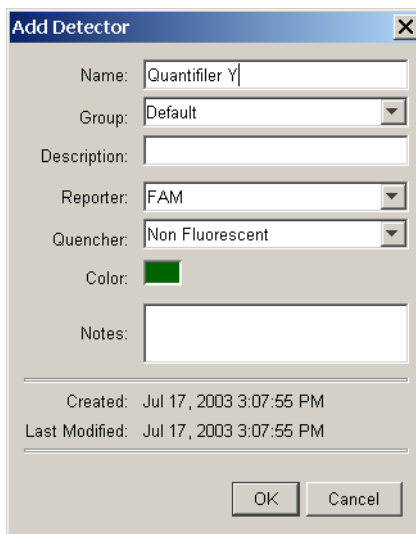

Add Detector


Name: Quantifiler Y

Group: Default

Description:

Reporter: FAM

Quencher: Non Fluorescent

Color: 

Notes:

Created: Jul 17, 2003 3:07:55 PM


Last Modified: Jul 17, 2003 3:07:55 PM

OK Cancel

- b. Click **OK** to return to the Detector Manager.

To create detectors: *(continued)*

4. Create a detector for the IPC assay:
- a. In the Detector Manager, click **New**, then complete the Add Detector dialog box:

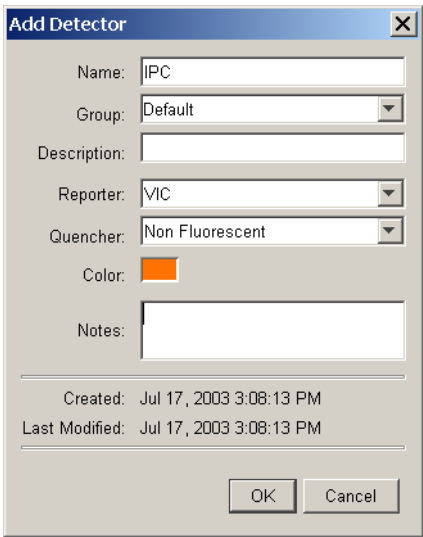

- b. Click **OK** to return to the Detector Manager.

**Copying
Detectors to the
Plate Document**

To copy detectors to the plate document:

1.	If the Detector Manager is not already open, select Tools > Detector Manager .
2.	Select the Quantifiler Human , Quantifiler Y , and the IPC detectors by clicking them while pressing the Ctrl key. Note: If the detectors are not available, create them first (see page 2-32 for the procedure).
3.	With the three detectors selected, click Copy To Plate Document .
4.	Click Done to close the Detector Manager and return to the plate window.

Applying Detectors for Standards

You need to apply the detectors to the plate document for the wells on the reaction plate that contain DNA quantification standards. Repeat the procedure until you complete applying detector tasks, quantities, and sample names for all quantification standards.

IMPORTANT! Set up detectors for each quantity and for each kit separately. For example, set up detectors for Std. 1 for the Quantifiler Human kit first, and then for Std. 2 for the Quantifiler Human kit, and so on, until you finish setting up the detectors for all wells containing quantification standards.

To apply detectors for quantification standards:

1.	In the plate grid, press the Ctrl key while you select the wells that correspond to a specific quantification standard for one kit.
2.	<p>Complete the Well Inspector:</p> <ol style="list-style-type: none"> Select the Use boxes for the applicable detectors: <ul style="list-style-type: none"> IPC Quantifiler Human or Quantifiler Y For the Quantifiler Human or Quantifiler Y detector: <ul style="list-style-type: none"> Click Unknown in the Task column, then select Standard from the drop-down list. Select the Quantity field and enter the quantity of DNA in the well. <p>IMPORTANT! Although you do not enter units for Quantity, you must use a consistent unit (for example, ng/μL) for all standard quantities. The units used for standard quantities defines the quantification units for analysis results.</p> <p>Note: Leave the IPC detector Task for standard reactions set to Unknown. Quantity values are not needed for IPC detectors.</p>

To apply detectors for quantification standards: *(continued)*

Step 2 continued:




- c. Enter the Sample Name (for example, Std. 1, Std. 2, and so on).
- d. Make sure that **ROX** is selected for the Passive Reference.

For example:

Setup | Instrument |

Well(s): A1-A2 ?

Sample Name: Std. 1

Use	Detector	Reporter	Task	Quantity	Color
<input checked="" type="checkbox"/>	IPC	VIC	Unknown	0	
<input checked="" type="checkbox"/>	Quantifiler Human	FAM	Standard	5E1	
<input type="checkbox"/>	Quantifiler Y	FAM		0	

Task for IPC set to **Unknown**
(default)

Applying
Detectors for
Unknown
Samples

You need to apply detectors to the plate document for the wells on the reaction plate that contain unknown samples.

IMPORTANT! If you run reactions for the Quantifiler Human kit and the Quantifiler Y kit on the same plate, apply detectors for unknown samples for each kit separately.

To apply detectors for unknown samples:




1. In the plate grid, press the **Ctrl** key and select the wells that contain unknown samples for one kit.
2. In the Well Inspector, select the Use boxes for the detectors in the selected wells:
 - IPC
 - Quantifiler Human or Quantifiler Y

For example:

Setup | Instrument |


Well(s): B5-D12 ?

Sample Name: * Mixed *

Use	Detector	Reporter	Task	Quantity	Color
<input checked="" type="checkbox"/>	IPC	VIC	Unknown	0	
<input checked="" type="checkbox"/>	Quantifiler Human	FAM	Unknown	0	
<input type="checkbox"/>	Quantifiler Y	FAM		0	

To apply detectors for unknown samples: *(continued)*

- | | |
|----|---|
| 3. | In the Well Inspector, make sure that ROX is selected for the Passive Reference. |
|----|---|

Passive Reference: ROX 

Adding Sample Names for Unknown Samples

Repeat this procedure to enter the names for all unknown samples.

To add sample names for unknown samples:

- | | |
|----|---|
| 1. | In the plate grid, select a reaction well containing an unknown sample. |
| 2. | <p>In the Well Inspector panel, enter a name in the Sample Name field.</p> <p>Note: Samples with identical sample names are treated as replicates by the SDS software. Results for replicate reactions are grouped together automatically for data analysis.</p> |


Setting Thermal
Cycler Conditions

Before running a Quantifiler kit assay, you need to make two changes to the default thermal cycler conditions:


- Thermal profile
- Sample volume

To set thermal cycler conditions:

1. In the plate window, select the **Instrument** tab.
2. Delete the Stage 1 hold step (50 °C for 2 minutes):
 - a. Press the **Shift** key and click within the Stage 1 hold step.

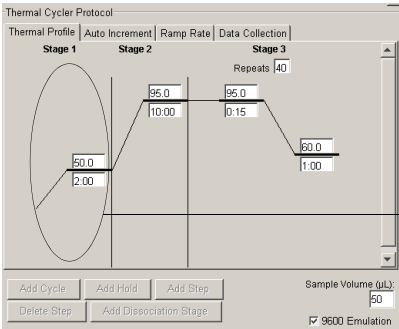


Press the Shift key and click within the Stage 1 hold step

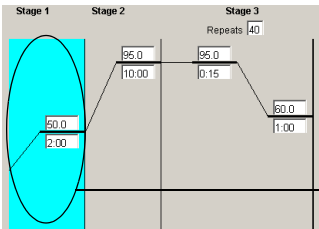


Hold step is selected

- b. After the hold step is selected, press the **Delete** key.

3. Make sure that the thermal profile appears as follows:

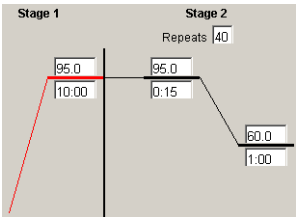

To set thermal cycler conditions: *(continued)*

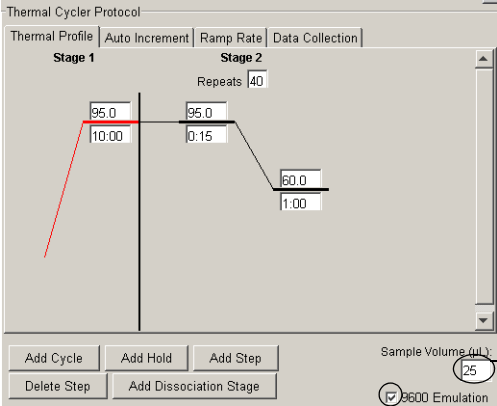4.	<p>Set the Sample Volume to 25 µL and make sure that the 9600 Emulation box is selected.</p> <p>Note: Selecting the 9600 Emulation box reduces the ramp rate.</p>  <p>Set the volume to 25 µL</p> <p>Make sure that this box is selected</p>
5.	<p>Make sure that the default settings are kept on the remaining tabs:</p> <ul style="list-style-type: none">• Auto Increment• Ramp Rate• Data Collection

Saving the Plate Document

Before running the reaction plate, save the plate document as an ABI PRISM SDS Single Plate (*.sds) file.

Note: To save the document as a template, see [“Setting Up a Plate Document Template” on page 2-40](#).

To save the plate document:

1.	Select File > Save As .
2.	For Files of Type, select ABI PRISM SDS Single Plate (*.sds) .

To save the plate document: *(continued)*

3.	Navigate to where you want to save the plate document file.
4.	In the File Name field, enter a name for the plate document.
5.	Click Save .

Setting Up a Plate Document Template

Purpose A plate document template reduces the time required to set up a plate document. This section describes how to create an SDS Template Document (*.sdt) set up for running Quantifiler kit assays.

Template Settings In addition to plate document settings (assay and container), templates can contain:

- Assay-specific detectors
- Well assignments for quantification standards, with detectors, tasks, and quantity
- Well assignments for unknown samples, with detectors and tasks
- Instrument settings: thermal cycler conditions and reaction volume settings.


Creating a Plate Document Template This procedure assumes that you have created the detectors for running reactions using the Quantifiler kits ([page 2-32](#)).

To create a plate document template:

1.	If the SDS software is not already started, select Start > Programs > Applied Biosystems > SDS 2.0 .
----	--

To create a plate document template: *(continued)*

2. Select **File > New**, then complete the New Document dialog box:

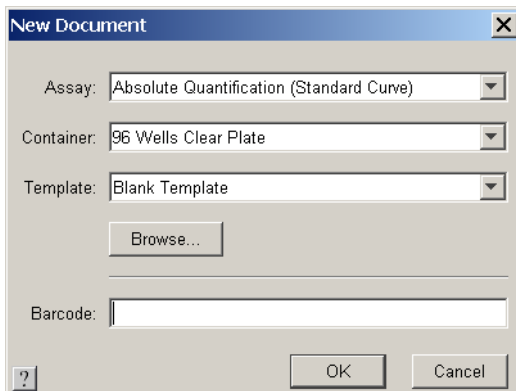

3. Apply the desired template settings to the plate document:
 - Copy detectors ([page 2-34](#))
 - Apply detectors for standards ([page 2-35](#))
 - Apply detectors for unknown samples ([page 2-36](#))
 - Set thermal cycler conditions ([page 2-38](#))
4. Select **File > Save As** and complete the Save As dialog box:
 - a. For Files of Type, select **ABI PRISM SDS Template Document (*.sdt)**.
 - b. Locate and select the Templates folder within the software folder:
X:Program Files > Applied Biosystems>7900HTSDS>Templates, where *X* is the hard drive on which the SDS software is installed.

Note: Saving the template file in the Templates folder makes it available in the Template drop-down list of the New Document dialog box (see [step 2](#) in “[Creating a Plate Document from a Template](#)” on [page 2-42](#)).
 - c. Enter a name for the template. For example, enter **Quantifiler Template**.
 - d. Click **Save**.

Creating a Plate Document from a Template

After you create a template, you can use it to create a plate document.

To create a plate document from a template:

1.	If the SDS software is not already started, select Start > Programs > Applied Biosystems > SDS 2.0 .
2.	<p>Select File > New and in the New Document dialog box and make the following selections:</p> <ul style="list-style-type: none"> • For Assay, select Absolute Quantitation. • For Container, select 96-Well Clear Plate. • For Template, select an appropriate template from the list. <p>Note: If the template is not available in the list, click Browse to locate and select an appropriate template.</p>
3.	<p>Complete the plate document setup:</p> <ul style="list-style-type: none"> • Copy detectors (page 2-34) • Apply detectors for standards (page 2-35) • Apply detectors for unknown samples (page 2-36) • Set thermal cycler conditions (page 2-38) <p>Note: The tasks that you perform vary according to which settings were defined in the template.</p>
4.	<p>Save the plate document (page 2-39).</p> <p>Note: For Files of Type, select ABI PRISM SDS Single Plate (*.sds).</p>

Chapter 3

PCR Amplification

PCR Amplification

This chapter covers:

Preparing the DNA Quantification Standard	3-2
Preparing the Reactions.	3-5
Running the Reactions	3-7

Preparing the DNA Quantification Standard

Required Materials

- Pipettors
- Pipette tips
- Quantifiler® Human DNA Standard

Note: The same standard can be used for both Quantifiler® kits.

- $T_{10}E_{0.1}$ buffer:
 - 10 mM Tris-HCl (pH 8.0)
 - 0.1 mM Na_2EDTA
 - 20 μ g/mL glycogen (optional)

Note: If you use $T_{10}E_{0.1}$ buffer with glycogen, you can store the DNA quantification standards for up to 2 weeks at 2 to 8 °C.

Guidelines for Calculating the Standards Dilution Series

The standard dilution series example shown in [Table 3-1 on page 3-3](#) is suitable for general use.

IMPORTANT! Applied Biosystems recommends:

- Three-fold dilution series with eight concentration points in the standard series for each assay
- Minimum input volume of 10 μ L DNA for dilutions (to ensure accuracy of pipetting)

Standards Dilution Series Example

Table 3-1 shows an example of one standards dilution series with the concentrations ranging from 50 ng/μL (Std. 1) to 0.023 ng/μL, or 23 pg/μL (Std. 8). A sample at the lowest concentration (2 μL per reaction) contains on average 14 to 16 copies of a diploid single-copy locus and 7 to 8 copies of a haploid single-copy locus.

Table 3-1 Standards dilution series example

Standard	Concentration (ng/μL)	Example Amounts	Minimum Amounts	Dilution Factor
Std. 1	50.000	50 μL [200 ng/μL stock] + 150 μL T ₁₀ E _{0.1} /glycogen buffer	10 μL [200 ng/μL stock] + 30 μL T ₁₀ E _{0.1} buffer	4×
Std. 2	16.700	50 μL [Std. 1] + 100 μL T ₁₀ E _{0.1} /glycogen buffer	10 μL [Std. 1] + 20 μL T ₁₀ E _{0.1} buffer	3×
Std. 3	5.560	50 μL [Std. 2] + 100 μL T ₁₀ E _{0.1} /glycogen buffer	10 μL [Std. 2] + 20 μL T ₁₀ E _{0.1} buffer	3×
Std. 4	1.850	50 μL [Std. 3] + 100 μL T ₁₀ E _{0.1} /glycogen buffer	10 μL [Std. 3] + 20 μL T ₁₀ E _{0.1} buffer	3×
Std. 5	0.620	50 μL [Std. 4] + 100 μL T ₁₀ E _{0.1} /glycogen buffer	10 μL [Std. 4] + 20 μL T ₁₀ E _{0.1} buffer	3×
Std. 6	0.210	50 μL [Std. 5] + 100 μL T ₁₀ E _{0.1} /glycogen buffer	10 μL [Std. 5] + 20 μL T ₁₀ E _{0.1} buffer	3×
Std. 7	0.068	50 μL [Std. 6] + 100 μL T ₁₀ E _{0.1} /glycogen buffer	10 μL [Std. 6] + 20 μL T ₁₀ E _{0.1} buffer	3×
Std. 8	0.023	50 μL [Std. 7] + 100 μL T ₁₀ E _{0.1} /glycogen buffer	10 μL [Std. 7] + 20 μL T ₁₀ E _{0.1} buffer	3×

Preparation Guidelines

While preparing the standards, keep in mind that:

- DNA quantification standards are critical for accurate analysis of run data
- Any mistakes or inaccuracies in making the dilutions directly affect the quality of results
- The quality of pipettors and tips and the care used in measuring and mixing dilutions affect accuracy

Preparing the DNA Quantification Standards

If you use $T_{10}E_{0.1}$ Buffer:

- With glycogen, you can store the prepared DNA quantification standards for up to 2 weeks at 2 to 8 °C.
- Without glycogen, long-term stability of the prepared DNA quantification standards may not be assured.

To prepare the DNA quantification standards dilution series:

1.	Label eight microcentrifuge tubes: Std. 1, Std. 2, Std. 3, and so on.
2.	Dispense the required amount of diluent ($T_{10}E_{0.1}$ Buffer with or without glycogen) to each tube.
3.	Prepare Std. 1: <ol style="list-style-type: none"> Vortex the Quantifiler Human DNA Standard 3 to 5 seconds. Using a new pipette tip, add the calculated amount of Quantifiler Human DNA Standard to the tube for Std. 1. Mix the dilution thoroughly.
4.	Prepare Std. 2 through 8: <ol style="list-style-type: none"> Using a new pipette tip, add the calculated amount of the prepared standard to the tube for the next standard. Mix the standard thoroughly. Repeat steps 4a and 4b until you complete the dilution series.

Preparing the Reactions

Required Materials

- Quantifiler® Human Primer Mix or Quantifiler® Y Human Male Primer Mix
- Quantifiler® PCR Reaction Mix
- 10-mL polypropylene tube
- 96-well reaction plate
- Extracted DNA samples
- DNA quantification standards dilutions series
- T₁₀E_{0.1} Buffer (with or without glycogen for negative controls)
- Optical Adhesive Cover

Preparing the Reactions

While preparing the reactions, keep the 96-well reaction plate in its base and do not place it on the counter.

To prepare the reactions:

1. Calculate the volume of each component needed to prepare the reactions, using the table below.

Component	Volume Per Reaction (μL)
Quantifiler Human Primer Mix or Quantifiler Y Human Male Primer Mix	10.5
Quantifiler PCR Reaction Mix	12.5

Note: Include additional reactions in your calculations to provide excess volume for the loss that occurs during reagent transfers.



CAUTION CHEMICAL HAZARD. Quantifiler PCR Reaction Mix may cause eye and skin irritation. Exposure may cause discomfort if swallowed or inhaled. Read the MSDS, and follow the handling instructions. Wear appropriate protective eyewear, clothing, and gloves.

To prepare the reactions: (continued)

2.	<p>Prepare the reagents:</p> <ul style="list-style-type: none">• Thaw the primer mix completely, then vortex 3 to 5 seconds and centrifuge briefly before opening the tube.• Swirl the Quantifiler PCR Reaction Mix gently before using. Do not vortex it.
3.	<p>Pipette the required volumes of components into an appropriately sized polypropylene tube.</p>
4.	<p>Vortex the PCR mix 3 to 5 seconds, then centrifuge briefly.</p>
5.	<p>Dispense 23 μL of the PCR mix into each reaction well.</p>
6.	<p>Add 2 μL of sample, standard, or control to the appropriate wells. For plate setup examples, see page 2-8, page 2-9, page 2-29, and page 2-30.</p> <p>IMPORTANT! Applied Biosystems recommends running duplicates of the eight DNA quantification standards for each assay and on each reaction plate (see page 3-4).</p>
7.	<p>Seal the reaction plate with the Optical Adhesive Cover.</p>
8.	<p>Centrifuge the plate at 3000 rpm for about 20 seconds in a tabletop centrifuge with plate holders to remove any bubbles.</p> <p>Note: If a tabletop centrifuge with 96-well plate adapters is not available, this step can be omitted.</p>
9.	<p>If you are using a 7000 or 7900HT instrument, place the compression pad over the Optical Adhesive Cover with the gray side down and the brown side up and with the holes positioned directly over the reaction wells.</p> <p>IMPORTANT! Do not use a compression pad if you are using a 7500 instrument.</p>

Running the Reactions

Before You Run the Reactions

Before you run the reactions, make sure that you have:

- Powered on the SDS instrument, computer, and software:
 - For 7000 SDS setup procedures, see [page 2-5](#)
 - For 7900HT SDS setup procedures, see [page 2-27](#)
- Set up a plate document for the run:
 - For 7000 SDS software procedures, see [page 2-7](#)
 - For 7900HT SDS software procedures, see [page 2-31](#)

Running the Plate on the 7000 SDS

To run the plate on the 7000 SDS:


1. Lift the handle at the bottom of the door on the front of the instrument until the door is raised completely. Gently push the carriage back until it stops and locks into place.



2. Position the plate in the instrument thermal block so that:
 - Well A1 is in the upper-left corner
 - The notched corner of the plate is in the upper-right corner

To run the plate on the 7000 SDS: (continued)

3. Gently push then release the carriage to unlatch it. The carriage automatically slides forward into position over the sample plate.

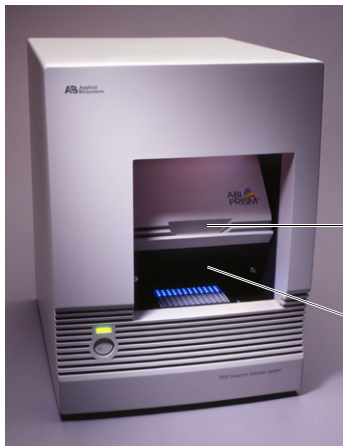

Do not pull the door forward by the handle

Gently push carriage back and release

4. After the door moves to the front, pull the handle down into place to close the cover.



CAUTION Do not pull the door handle to move the carriage forward. This may cause serious damage to the door or the door mechanism.



WARNING PHYSICAL INJURY HAZARD.

During instrument operation, the temperature of the heated cover can be as high as 108 °C, and the temperature of the sample block can be as high as 100 °C. Before performing the procedure, keep hands away until the heated cover and sample block reach room temperature.

5. In the SDS software, open the plate document that you set up for the run.
6. Select the **Instrument** tab, then click **Start**.

Running the Plate on the 7900HT SDS

To run the plate on the 7900HT SDS:

1.	In the SDS software, select the Instrument tab for the plate document.
2.	In the Real-Time tab, click Open/Close to rotate the instrument tray to the OUT position.
3.	Place the plate in the instrument tray so that: <ul style="list-style-type: none"> Well A1 is in the upper-left corner The notched corner is in the upper-right corner
4.	Click Start to rotate the instrument tray to the IN position and to start the run. Note: The instrument may pause to allow the heated cover to heat to the appropriate temperature before beginning the run. The SDS software collects and saves the run data and the Real-Time tab displays the instrument status and run progress.
5.	After the run is complete, remove the plate from the instrument: <ol style="list-style-type: none"> Click Open/Close in the Instrument tab of the plate document that is open and connected to the 7900HT instrument. The instrument tray rotates to the OUT position. Remove the plate from the instrument. Click Open/Close in the Instrument tab to rotate the instrument tray to the IN position.

Chapter 4

Data Analysis and Results

Data Analysis and Results

This chapter covers:

Section 4.1 7000 SDS Data Analysis	4-3
Analyzing the Plate Document	4-3
Viewing Results.	4-4
Section 4.2 7900HT SDS Data Analysis.	4-7
Analyzing the Plate Document	4-7
Viewing Results.	4-8


Section 4.1 7000 SDS Data Analysis

Analyzing the Plate Document

Analyze a run after it is complete and reanalyze after you make any changes to the plate document, such as sample names.

To analyze a plate document:

1. Open the plate document to analyze.
2. Verify the analysis settings:
 - a. Select **Analysis > Analysis Settings** to open the Analysis Settings dialog box.
 - b. Verify that the settings are as shown below, then click **OK**.

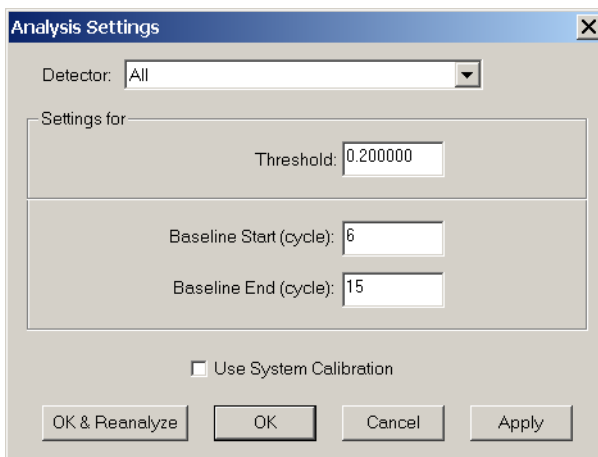

IMPORTANT! If the analysis settings differ from those shown here, change them to match the settings before clicking OK.

3. Select **Analysis > Analyze**.

Viewing Results

Overview Viewing the results of data analysis can involve one or more of the following:

- [Viewing the Standard Curve \(page 4-4\)](#)
- [Viewing the Amplification Plot \(page 4-5\)](#)
- [Viewing the Report \(page 4-5\)](#)
- [Printing or Exporting the Report \(page 4-6\)](#)

Viewing the Standard Curve

For information about interpreting and troubleshooting the standard curve, see [“Examining the Standard Curve” on page 5-4](#) and [“Troubleshooting the Standard Curve” on page 5-6](#).

To view the standard curve:

1.	In the Results tab, select the Standard Curve tab.
2.	In the Detector drop-down list, select the detector that corresponds to the kit that you are using: <ul style="list-style-type: none">• Quantifiler Human <i>or</i>• Quantifiler Y
3.	View the C_T values for the quantification standard reactions and the calculated regression line, slope, y-intercept, and R^2 values.

Amplification Plot Results

The amplification plot can display one of the following:

- Plot of normalized reporter signal (R_n) versus cycle number for each reaction
- C_T versus well position on the assay plate

For more information about the amplification plot, see [“Real-Time Data Analysis” on page 1-10](#).

Viewing the Amplification Plot

For troubleshooting information, see [“Troubleshooting Amplification Plots” on page 5-12.](#)

To view the amplification plot:

1.	In the Results tab, select the Amplification Plot tab.
2.	In the Detector drop-down list, select a detector: <ul style="list-style-type: none"> • Quantifiler Human or Quantifiler Y • IPC
3.	Select the appropriate samples in the table below the amplification plot.
4.	Make sure that the Threshold is set to 0.20 , the default setting. Note: If you move the threshold bar, it changes from green to red to indicate reanalysis is needed. After reanalysis, it changes from red to green.

Viewing the Report

The report summarizes the quantity of DNA present in the samples. For information about the quantities reported, see [“Assessing Quantity” on page 5-16.](#)

To view the report:


1.	In the analyzed plate document, select the Results tab, then select the Report tab.
2.	Select the reactions in the 96-well plate representation below the report to display the results in the report.
3.	View the Qty column to determine the quantity of DNA in each sample. Note: Quantities are calculated only if quantification standards were run and set up correctly in the software. Otherwise, only C_T values are shown.

Printing or
Exporting the
Report

For more information about exporting data, see the *ABI PRISM® 7000 Sequence Detection System User Guide* (PN 4330228).

To print or export the report:

1. In the Report tab of the Results window, select **Tools > Report Settings**, then set up how the report is printed and exported:

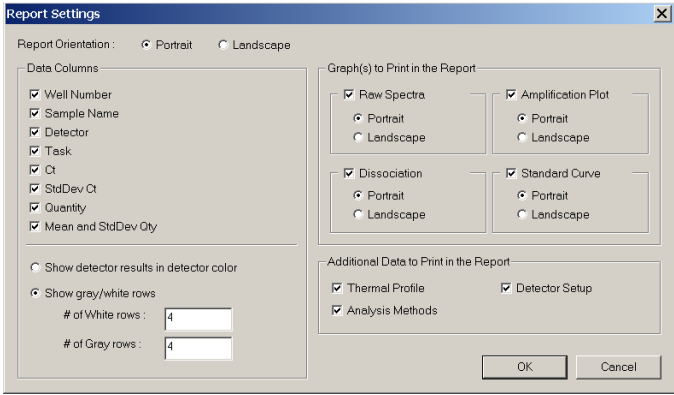

- 2.
- Select **File > Print** to print the report.
 - Select **File > Export** to export the report as tab-delimited text.

Note: You can later open the exported file using spreadsheet software.


Section 4.2 7900HT SDS Data Analysis

Analyzing the Plate Document

Analyze a run after it is complete and reanalyze the run:

- Each time that you open a plate document to convert the saved raw data into analyzed data
- After you make changes to the plate document, such as sample names

To analyze the plate document:

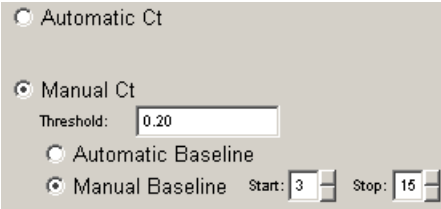1.	Open the plate document to analyze.
2.	<p>Verify the analysis settings:</p> <ol style="list-style-type: none"> Select Analysis > Analysis Settings to open the Analysis Settings dialog box. Verify that the settings are as shown below, then click OK:  <p>IMPORTANT! If the analysis settings differ from those shown here, change them to match the settings before clicking OK.</p>
3.	Select Analysis > Analyze for the software to convert the raw data to analyzed data.
4.	Select the Results tab to view the results.

Viewing Results

Overview Viewing the results of data analysis can involve one or more of the following:

- [Viewing the Standard Curve \(page 4-8\)](#)
- [Viewing the Amplification Plot \(page 4-9\)](#)
- [Viewing the Results Table \(page 4-9\)](#)
- [Printing the Results \(page 4-10\)](#)
- [Exporting the Results \(page 4-10\)](#)

Viewing the Standard Curve For information about interpreting and troubleshooting the standard curve, see [“Examining the Standard Curve” on page 5-4](#) and [“Troubleshooting the Standard Curve” on page 5-6](#).

To view the standard curve:

1.	In the Results tab, select the Standard Curve tab.
2.	In the Detector drop-down list, select the detector that corresponds to the kit that you are using: <ul style="list-style-type: none">• Quantifiler Human <i>or</i>• Quantifiler Y
3.	View the C_T values for the quantification standard reactions and the calculated regression line, slope, intercept, and R^2 values.

Amplification Plot Results

The amplification plot can display one of the following:

- Plot of normalized reporter signal (R_n) versus cycle number for each reaction
- C_T versus well position on the assay plate

For more information about the amplification plot, see [“Real-Time Data Analysis” on page 1-10](#).

Viewing the Amplification Plot

For troubleshooting information, see [“Troubleshooting Amplification Plots” on page 5-12.](#)

To view the amplification plot:

1.	After the run is finished, select the Results tab, then select the Amplification Plot tab.
2.	In the Detector drop-down list, select the detector: <ul style="list-style-type: none"> • Quantifiler Human or Quantifiler Y • IPC
3.	Select the appropriate samples in the 96-well grid or the sample table to the left of the amplification plot.
4.	Make sure that the Threshold is set to 0.20 , the default setting.

Results Table

The results table displays:

- Well position of samples
- Sample names
- Detector assignments
- Task assignments
- C_T values
- Quantity
- Mean and standard deviation for C_T values and Quantity, if replicate groups were defined in assay setup

Viewing the Results Table

View the Qty column to determine the quantity of DNA present in each sample.

Note: Units for calculated quantities are not displayed but are the same as those specified for the quantification standards when you set up the plate document.

Note: Quantities are calculated only if quantification standards were run and set up correctly in the software. Otherwise, only C_T values are shown.

For more information about the quantities reported, see [“Assessing Quantity” on page 5-16.](#)

Printing the Results

To print the results:

1.	Select File > Print Report .
2.	Select the data to include in the report by selecting the corresponding boxes for: <ul style="list-style-type: none"> • Document Information • Thermal Cycler Conditions • Detector Information • Well Status Summary • Raw Data Plot • Multicomponent Data Plot • Amplification Plot
3.	Click Page Setup , then select: <ul style="list-style-type: none"> • Header/footer information and placement • Layout orientation and size
4.	Click Print to print the report.

Exporting the Results

You can export the results in tab-delimited (*.txt) format and later open the exported files using spreadsheet software.

To export the results:

1.	Select File > Export .
2.	Select the results to export: <ul style="list-style-type: none"> • Setup Table • Results Table • Multicomponent • Clipped
3.	Select whether you want to export data from all wells or selected wells.
4.	Select the SDS format of data to export.
5.	Select Group by replicates if you want the replicates to be grouped together in the exported results.

To export the results: (continued)

6.	Locate, then select the folder where you want to save the exported results file.
7.	Enter the File Name, then click Export .

Interpretation of Results

This chapter covers:

Checking Analysis Settings	5-2
Examining the Standard Curve	5-4
Troubleshooting the Standard Curve	5-6
Using the Internal PCR Control System	5-10
Troubleshooting Amplification Plots	5-12
Assessing Quantity	5-16

Checking Analysis Settings


The validity of the results requires correct analysis settings.

Checking Analysis Settings on the 7000 SDS

To check analysis settings on the 7000 SDS:

1. If the SDS software is not already started, select **Start > Programs > ABI Prism 7000 > ABI Prism 7000 SDS Software**.
2. Select **File > Open**.
3. Locate the plate document for the assay run of interest, select it, then click **Open**.
4. Select **Analysis > Analysis Settings**.


5. For all detectors, confirm that the settings are as shown below:

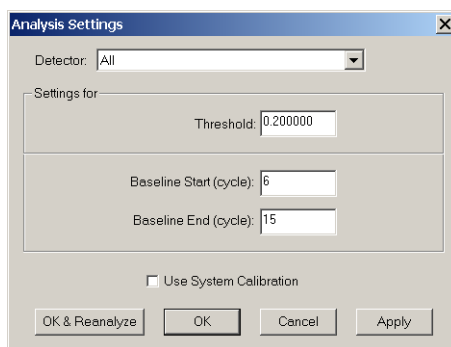

6. If the analysis settings differ from those shown in [step 5](#):
 - a. Change the settings to match those in [step 5](#).
 - b. Click **Apply**.
 - c. Click **OK & Reanalyze** to close the dialog box and reanalyze the plate document.
 - d. View the results using [Chapter 4, “Data Analysis and Results.”](#)

Checking Analysis Settings on the 7900HT SDS

To check analysis settings on the 7900HT SDS:

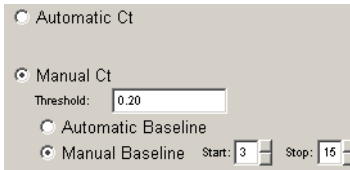1.	If the SDS software is not already started, select Start > Programs > Applied Biosystems > SDS 2.0 .
2.	Select File > Open .
3.	Locate the plate document for the assay run of interest, select it, then click Open .
4.	<p>Select Analysis > Analysis Settings and confirm that the settings are as shown below:</p> 
5.	<p>If the analysis settings differ from those shown in step 4:</p> <ol style="list-style-type: none"> Change the settings to match those in step 4. Click OK. Select Analysis > Analyze for the software to reanalyze the data. View the results using Chapter 4, “Data Analysis and Results.”

Examining the Standard Curve

Examine the standard curve results to evaluate the quality of the results from the quantification standard reactions.

About Standard Curve Results

The standard curve is a graph of the C_T of quantification standard reactions plotted against the starting quantity of the standards. The software calculates the regression line by calculating the best fit with the quantification standard data points. The regression line formula has the form:

$$C_T = m [\log (Qty)] + b$$

where m is the slope, b is the y-intercept, and Qty is the starting DNA quantity. The values associated with the regression analysis can be interpreted as follows:

- **R² value** – Measure of the closeness of fit between the standard curve regression line and the individual C_T data points of quantification standard reactions. A value of 1.00 indicates a perfect fit between the regression line and the data points.
- Regression coefficients:
 - **Slope** – Indicates the PCR amplification efficiency for the assay. A slope of -3.3 indicates 100% amplification efficiency.
 - **Y-intercept** – Indicates the expected C_T value for a sample with Qty = 1 (for example, 1 ng/ μ L)

R² Value

An R^2 value ≥ 0.99 indicates a close fit between the standard curve regression line and the individual C_T data points of quantification standard reactions

If the R^2 value is <0.98 check the following:

- Quantity values entered for quantification standards in the Well Inspector during plate document setup
- Making of serial dilutions of quantification standards
- Loading of reactions for quantification standards
- Failure of reactions containing quantification standards
- C_T value for Std. 8 of the DNA quantification standard (23 pg/ μ L), if using the Quantifiler[®] Y kit

R² Value < 0.98 for Quantifiler Y Kit Only

If the R² value is <0.98 for the Quantifiler Y kit only, you may choose to omit Std. 8 of the DNA quantification standard (23 pg/μL) from analysis.

At the lowest concentration point, there are only 7 to 8 copies per 2 μL reaction of the haploid target locus for the Quantifiler Y kit. Because of stochastic effects when using the lowest concentration point with Quantifiler Y kits, the C_T values are more variable at the lowest concentration point and may affect the closeness of fit between the standard curve regression line and the individual data points of the quantification standard.

To omit Std. 8 from analysis (for Quantifiler Y kits only):

1.	Select the wells in the plate document that correspond to Std. 8 and open the Well Inspector.
2.	Change the Task assignment for the Quantifiler Y detector from Standard to Unknown .
3.	Reanalyze the plate to incorporate the change.

Slope

A slope close to –3.3 indicates optimal, 100% PCR amplification efficiency.

Table 5-1 Range and average of standard curve slope values

Kit	Typical Slope (range)	Average Slope
Quantifiler Human	–2.9 to –3.3	–3.1
Quantifiler Y	–3.0 to –3.6	–3.3

If the slope varies beyond the typical range indicated in [Table 5-1](#), check the following:

- Assay setup
- Software setup
- Reagents
- Instrument

Troubleshooting the Standard Curve

Table 5-2 Troubleshooting the standard curve

Observation	Possible Cause	Recommended Action
<p>Slope for the standard curve differs significantly from -3.33 or</p> <p>R^2 value significantly less than 0.98 to 0.99</p>	When applying detectors for standards, the Task and Quantity were applied to the wrong detector (see “Example 1” on page 5-7).	<ol style="list-style-type: none"> 1. From the plate document, double-click a well containing a DNA quantification standard to view the Well Inspector. 2. Verify that the Task and Quantity were applied to the correct detector and reanalyze.
	When applying detectors for the standards, the incorrect Quantity was entered (see “Example 2” on page 5-8).	<ol style="list-style-type: none"> 1. From the plate document, double-click a well containing a DNA quantification standard to view the Well Inspector. 2. Verify that the correct Quantity was entered and reanalyze.
	Stochastic effects when using the lowest concentration point with the Quantifiler Y kit.	Omit Std. 8 of the DNA quantification standard (23 pg/ μ L) from analysis.
<p>At each concentration in the standard curve:</p> <ul style="list-style-type: none"> • There are four replicates • There is a large difference in C_T between the replicates <p>Note: This observation applies only when Quantifiler Human kit reactions and Quantifiler Y kit reactions are run together on the same reaction plate.</p>	The same detector was applied for the Quantifiler Human kit standard reactions and for the Quantifiler Y kit standard reactions (see “Example 3” on page 5-9).	<ol style="list-style-type: none"> 1. From the plate document, double-click a well containing a DNA quantification standard to view the Well Inspector. 2. Verify that the correct detector is in use and that the Task and Quantity were applied to the correct detector and reanalyze.


The examples shown in the following sections can be caused by errors made in applying the detectors for standards when setting up the plate document. For instructions on how to apply the detectors for standards, see:

- Page 2-15 (7000 SDS)
- Page 2-35 (7900HT SDS)

Note: The standard curves shown in these examples are not optimal and should not be used.


Example 1 Observation

Almost all of the C_T values for the DNA quantification standard reactions lie outside of the standard curve and form a straight horizontal line.

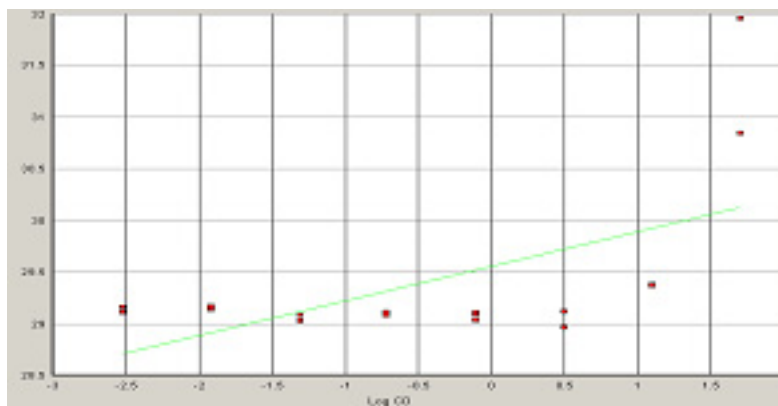


Possible Cause

When applying detectors for the standards, the Task and Quantity were applied to the IPC detector instead of to the Quantifiler Human detector.

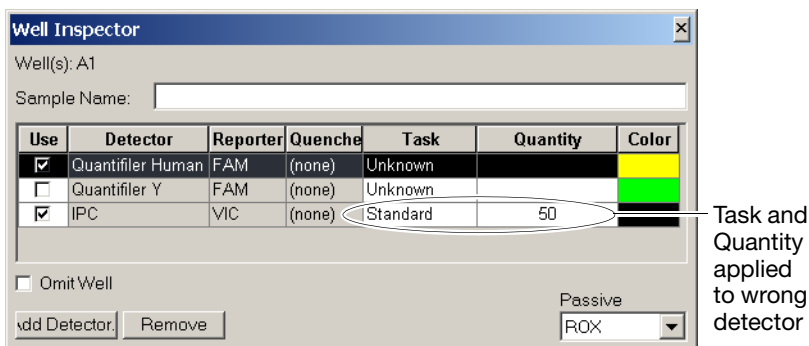

Example 2 Observation

One point lies outside of the standard curve.

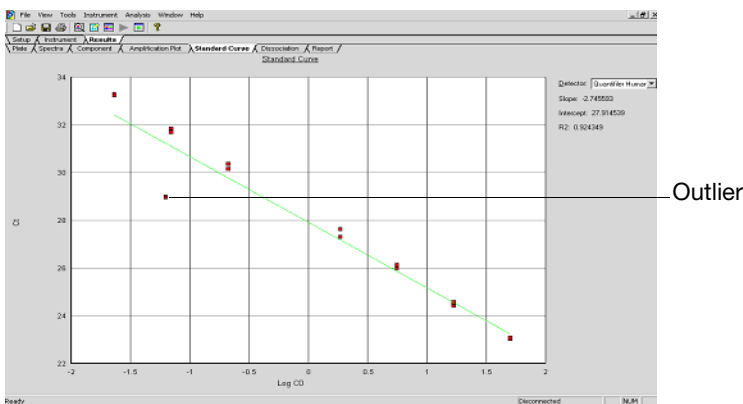

Possible Cause

When applying detectors for the standards, the incorrect Quantity was entered. In the example shown above, 0.062 was entered for the Quantity instead of 0.62.

Well Inspector

Well(s): E1

Sample Name:

Use	Detector	Reporter	Quencher	Task	Quantity	Color
<input checked="" type="checkbox"/>	Quantifiler Human	FAM	(none)	Standard	0.062	Yellow
<input type="checkbox"/>	Quantifiler Y	FAM	(none)	Unknown		Green
<input checked="" type="checkbox"/>	IPC	VIC	(none)	Unknown		Black

☐ Omit Well


Passive
 ROX

Incorrect Quantity entered for standard

Example 3 Observation

At each concentration in the standard curve:

- There are four replicates
- There is a large difference in the C_T between the replicates

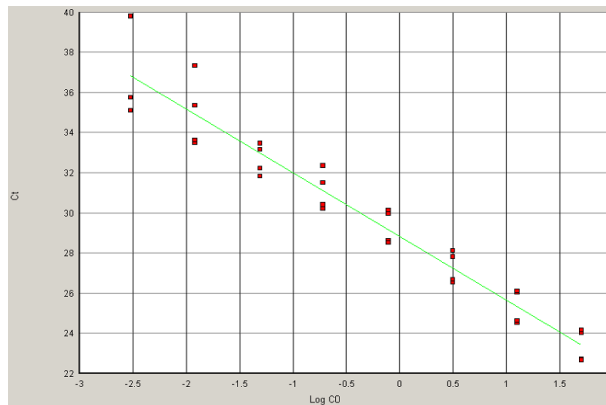

Possible Cause

The Quantifiler® Human kit assay and the Quantifiler® Y kit assay were performed on the same reaction plate and when applying detectors for standards, the same detector was applied for Quantifiler Human kit standard reactions and for the Quantifiler Y kit standard reactions.

Using the Internal PCR Control System

Purpose Use the Internal PCR Control (IPC) system to distinguish between true negative sample results and reactions affected by:

- The presence of PCR inhibitors
- Assay setup
- A chemistry or instrument failure

Components The following components of the IPC system are present in the Quantifiler PCR mix:

- Synthetic DNA template
- Primers that hybridize specifically to the synthetic DNA template
- Probe labeled with VIC[®] dye

Interpreting IPC Results In the amplification plot window of the SDS software, observe amplification of the FAM[™] dye (Quantifiler Human detector or Quantifiler Y detector) and the VIC[®] dye (IPC detector), then use [Table 5-3](#) to interpret the IPC results.

Table 5-3 Interpreting IPC amplification results

Quantifiler Human or Quantifiler Y (FAM Dye)	IPC (VIC Dye)	Interpretation
No amplification	Amplification	True negative
No amplification	No amplification	Invalid result
Amplification (low C_T and high ΔR_n)	No amplification	Disregard IPC result
Amplification (high C_T and low ΔR_n)	No amplification	Partial PCR inhibition

Note: Positive amplification is when the C_T value for the detector is <40 . Because samples contain unknown amounts of DNA, a large range of C_T values is possible. Because the IPC system template DNA is added to the reaction at a fixed concentration, the $C_{T\text{ VIC}}$ should range from 20 to 30.

True Negative Results	<p>With a true negative result:</p> <ul style="list-style-type: none">• FAM dye signal indicates that the human-specific target failed to amplify• VIC dye signal ($C_{T\text{ VIC}}$ between 20 and 30) indicates that the IPC target was amplified; so, the PCR was not inhibited
Invalid IPC Results	<p>If the human-specific target <i>and</i> the IPC target failed to amplify, it is not possible to distinguish between the absence of DNA and PCR inhibition.</p>
Disregarding IPC Results	<p>With extremely high concentrations of human genomic DNA (>10 ng/μL), competition between the human-specific and IPC PCR reactions appears to suppress IPC amplification for that sample. If the target amplifies with low C_T and high ΔR_n results, it is unlikely that PCR inhibitors are present. In these cases, appearance of suppression or failure of IPC amplification can be disregarded.</p>
Partial PCR Inhibition	<p>Weak amplification (high C_T value and low ΔR_n value) of the human target and no amplification of the IPC may indicate partial PCR inhibition in the sample.</p>
Determining the Normal Range for IPC	<p>To determine the normal range of C_T values for the IPC, view the VIC dye signal in the amplification plots for the quantification standards. If the assays were set up properly and the buffer used to dilute the quantification standards was free of PCR inhibitors, the reactions should show normal IPC amplification across a broad range of input DNA.</p>
Evaluating PCR Inhibition	<p>If the IPC amplification for certain samples appears reduced relative to IPC amplification for quantification standards, the decreased IPC amplification may be interpreted as partial PCR inhibition. The IPC results can help you decide the next step:</p> <ul style="list-style-type: none">• Proceed directly to an STR assay of the sample• Repeat the DNA extraction from the sample• Perform additional cleanup of the sample

Troubleshooting Amplification Plots

Table 5-4 Troubleshooting amplification plots

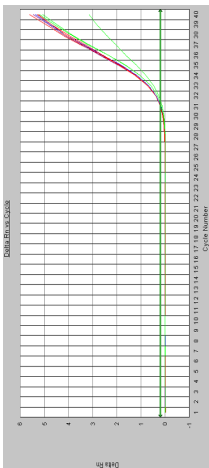
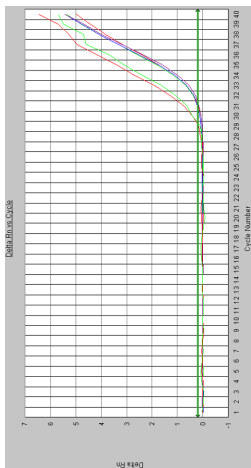
Observation	Possible Cause	Recommended Action
 <p>ΔR_n and C_T values inconsistent with replicates</p>	<p>Evaporation of reaction mixture from some wells because the Optical Adhesive Cover was not sealed to the reaction plate properly or the compression pad was not used during the run</p>	<p>Confirm the cause:</p> <ol style="list-style-type: none"> 1. Select the Component tab. Affected wells should generate significantly less fluorescence compared to unaffected replicates. 2. Check the amount of solution in each well of the reaction plate. Wells affected by evaporation should contain less solution compared to unaffected wells and should correspond with the inconsistent results. <p>For subsequent runs, make sure that the Optical Adhesive Cover is sealed to the reaction plate properly and that the compression pad is used.</p>
 <p>ΔR_n and C_T values inconsistent with replicates</p>	<p>Incorrect volume of Quantifiler PCR Reaction Mix added to some reactions</p>	<p>Confirm the cause:</p> <ol style="list-style-type: none"> 1. Select the Component tab. Affected wells should generate significantly different amounts of fluorescence compared to unaffected replicates. 2. Select the Spectra tab. Wells with the incorrect volume of Quantifiler PCR Reaction Mix should generate significantly different amounts of fluorescence compared to unaffected wells.

Table 5-4 Troubleshooting amplification plots (*continued*)


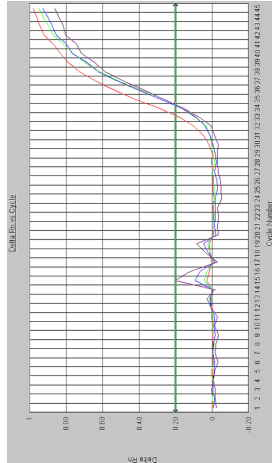
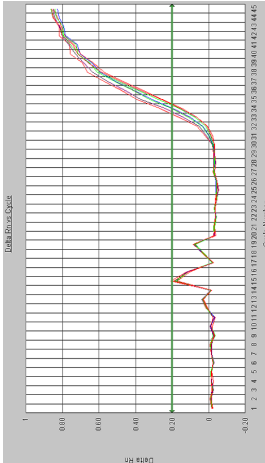
Observation	Possible Cause	Recommended Action
 <p>Jagged amplification plots</p>	Weak lamp or improper replacement	Replace the lamp or make sure that the lamp was replaced properly.
	Mechanical or optical misalignment	<ol style="list-style-type: none"> 1. Localize the wells that contain baseline spikes. 2. Run the TaqMan® RNase P Instrument Verification Plate (PN 4310982). 3. Perform the instrument function tests. <ul style="list-style-type: none"> – If a function test fails, contact your Applied Biosystems Service Representative. – If all functional tests pass, the reaction plate or the door of the instrument may not have been aligned properly during the run. <p>Note: See your instrument user guide for instructions on how to perform instrument function tests.</p>
 <p>Baseline spikes with certain reactions and normal amplification with other reactions</p>		

Table 5-4 Troubleshooting amplification plots (continued)

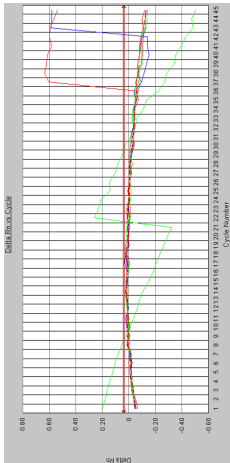
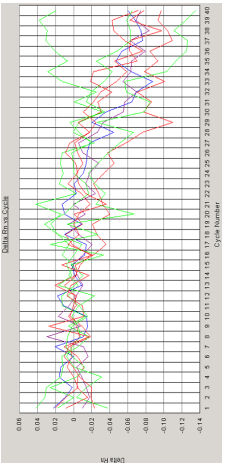
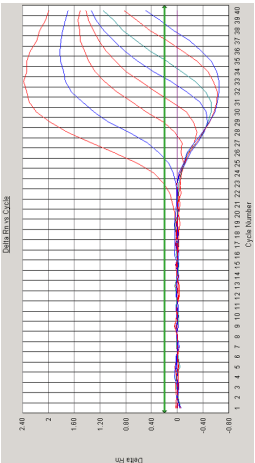
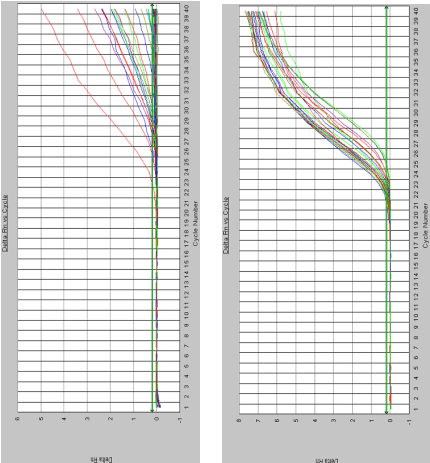
Observation	Possible Cause	Recommended Action
 <p>Abnormal amplification plots in one column of reactions</p>	<p>Uncalibrated pure dyes, damage to the lens, or dust on the filters and/or mirror</p>	<p>If the pure dyes are not calibrated, run the pure dyes and recalibrate.</p> <p>Note: See your instrument user guide for instructions on how to run pure dyes and recalibrate.</p>
 <p>No defined amplification plots</p>	<p>Incorrect detector selected on the amplification plot or incorrect detector applied to the reactions when setting up the plate document</p>	<ol style="list-style-type: none"> 1. Make sure that the correct detector is selected on the amplification plot. 2. If the amplification plots are still not defined: <ol style="list-style-type: none"> a. From the plate document, double-click a well to view the Well Inspector. b. Verify that the detector settings are correct and reanalyze.
 <p>Abnormal ΔR_n values and some negative R_n values</p>	<p>Incorrect passive reference was selected when setting up the plate document</p>	<p>Confirm the diagnosis:</p> <ol style="list-style-type: none"> 1. From the plate document, double-click a well to view the Well Inspector. 2. Observe which Passive Reference is selected. <p>Note: ROX should be selected as the Passive Reference.</p>

Table 5-4 Troubleshooting amplification plots (*continued*)

Observation	Possible Cause	Recommended Action
 <p>Reactions in rows B, C, and D show poor amplification and reactions in the rows E, F, and G show good amplification</p>	<p>Instrument door was not aligned properly on the reaction plate</p>	<ol style="list-style-type: none"> 1. Localize the wells that show poor amplification. 2. Run the TaqMan® RNase P Instrument Verification Plate (PN 4310982). 3. Check the calibration of the regions of interest (ROI). 4. Perform the instrument function tests. <ul style="list-style-type: none"> – If a function test fails, contact your Applied Biosystems Service Representative. – If all functional tests pass, the reaction plate or the door of the instrument may not have been aligned properly during the run. <p>Note: See your instrument user guide for instructions on how to check ROI calibration and to perform instrument function tests.</p>

Assessing Quantity

Purpose After viewing the results and assessing the quality of the results, the analyst should determine whether sufficient DNA is present to proceed with a short tandem repeat (STR) assay.

Assay Sensitivity Quantifiler kit assays can detect < 23 pg/ μ L of human genomic DNA in samples. For samples loaded at 2.0 μ L per reaction, this concentration corresponds to < 13 copies of the Quantifiler Human target DNA and < 7 copies of the Quantifiler Y target locus (Y chromosome loci are haploid).

Stochastic Effects In the 23-pg/ μ L concentration range, stochastic effects, or the statistical effect of sampling low-copy loci, may cause significant variability in assay results.

Validity The detection and quantification of low-copy DNA samples with the Quantifiler kits is valid. However, the amounts present in the sample may be below the working range of certain genotyping methods.

If Insufficient DNA Is Present If the results from Quantifiler kit reactions indicate that insufficient DNA is present to perform an STR assay, the analyst may decide to:

- Extract the DNA again, then repeat the test with the Quantifiler kit before performing STR analysis
- Concentrate the sample, then repeat the test with the Quantifiler kit before performing STR analysis

Experiments and Results

This chapter covers:

Overview	6-2
Section 6.1 ABI PRISM® 7000 Sequence Detection System Validation (SDS Software v1.0)	6-3
Section 6.2 Applied Biosystems 7900HT Real-Time PCR System Validation (SDS Software v2.0)	6-37
Section 6.3 Casework Sample Analysis	6-46
Section 6.4 Applied Biosystems 7500 Real-Time PCR System Validation (SDS Software v1.2.3)	6-50
Section 6.5 ABI PRISM® 7000 Sequence Detection System Validation (SDS Software v1.2.3)	6-65

Note: The information in Sections 6.4 and 6.5 is also contained in the *Quantifiler® Kits User Bulletin: Validation Using SDS Software Version 1.2.3 on the Applied Biosystems 7500 Real-Time PCR System and the ABI PRISM® 7000 Sequence Detection System* (PN 4374659 Rev. A, 4/2006).

Overview

About This Chapter

This chapter provides results of the validation experiments performed by Applied Biosystems using the Quantifiler® Human DNA Quantification Kit and the Quantifiler® Y Human Male DNA Quantification Kit.

Importance of Validation

Although the Quantifiler kits are not DNA genotyping assays, they are intended for use before performing genotyping assays such as the AmpFSTR® PCR Amplification kits. By testing the procedure with samples commonly encountered in forensic and parentage laboratories, the validation process clarifies attributes and limitations that are critical for sound data interpretation in casework.

Experiments

Experiments to evaluate the performance of the Quantifiler kits were performed at Applied Biosystems, according to the DNA Advisory Board (DAB) Quality Assurance Standards For Forensic DNA Testing Laboratories (DAB, 1998). These DAB standards describe the quality assurance requirements that a laboratory should follow to ensure the quality and integrity of the data and competency of the laboratory. The DAB defines a laboratory as a facility in which forensic DNA testing is performed. Additional validation was performed according to guidelines from the Scientific Working Group on DNA Analysis Methods (SWGDAM).

The experiments focused on kit performance parameters relevant to the intended use of the kits as human-specific DNA quantification assays and as a part of a forensic DNA genotyping procedure.

Each laboratory using the Quantifiler® Human DNA Quantification Kit or the Quantifiler® Y Human Male DNA Quantification Kit should perform appropriate validation studies.

Section 6.1 ABI PRISM® 7000 Sequence Detection System Validation (SDS Software v1.0)

DAB (DNA Advisory Board) Guideline 8.1.1

“Developmental validation that is conducted shall be appropriately documented.” (DAB, 1998).

DAB Guideline 8.1.2

“Novel forensic DNA methodologies shall undergo developmental validation to ensure the accuracy, precision and reproducibility of the procedure.” (DAB, 1998).

6.1.1 Precision6-4

6.1.2 Reproducibility6-7

DAB Guideline 8.1.2.2

“Species specificity, sensitivity, stability and mixture studies are conducted.” (DAB, 1998).

6.1.3 Specificity with a Human DNA Panel6-10

6.1.4 Specificity with a Non-Human Panel6-11

6.1.5 Specificity with a Bacterial Pools Panel6-14

6.1.6 Sensitivity6-16

6.1.7 Stability6-17

6.1.8 Mixture Studies6-21

6.1.9 Degraded DNA Studies6-23

6.1.10 Comparisons with Other Methods6-27

6.1.11 Comparison with A260 and Quantiblot Kit6-27

6.1.12 Comparison with A260 and Dye Intercalation6-31

6.1.13 Assay Background6-35

6.1.1 Precision

The precision of the Quantifiler® Human kit and the Quantifiler® Y kit was tested by performing runs on different instruments and on different days.

Experiment One set of eight serial dilutions of the Quantifiler Human DNA Standard was prepared. The dilutions ranged from 50 ng/μL to 23 pg/μL in three-fold increments.

Three different reaction plates were prepared and each plate contained duplicate reactions of the dilutions using the Quantifiler Human kit and the Quantifiler Y kit.

The three plates were run on three different 7000 SDS instruments, using standard thermal cycler conditions for the Quantifiler kits. The multiple runs were performed on two different days, using the same three 7000 SDS instruments.

The $C_{T\text{ FAM}}$ values were recorded and the means and standard deviations of the $C_{T\text{ FAM}}$ values were calculated for each of the eight dilutions using the Quantifiler Human kit and the Quantifiler Y kit.

Results Table 6-1 shows the means and standard deviations of the $C_{T\text{ FAM}}$ values calculated for all 12 reactions of each quantification standard dilution for the Quantifiler Human and Quantifiler Y kits.

Table 6-1 Precision: C_T values

Quantification Standard Dilution (ng/μL)	Quantifiler Human Kit		Quantifiler Y Kit	
	C_T (Mean)	Standard Deviation	C_T (Mean)	Standard Deviation
50	23.09	0.10	23.94	0.21
16.7	24.64	0.17	25.38	0.17
5.56	26.19	0.16	26.91	0.13
1.85	27.67	0.17	28.35	0.15
0.62	29.09	0.17	29.84	0.26
0.21	30.31	0.19	31.38	0.31

Table 6-1 Precision: C_T values (continued)

Quantification Standard Dilution (ng/ μ L)	Quantifiler Human Kit		Quantifiler Y Kit	
	C_T (Mean)	Standard Deviation	C_T (Mean)	Standard Deviation
0.068	31.90	0.28	33.38	0.44
0.023	33.45	0.48	35.19	0.73

Figure 6-1 and Figure 6-2 show the $C_{T\text{ FAM}}$ results for all 8 quantification standard dilutions reactions using the Quantifiler Human kit and the Quantifiler Y kit.

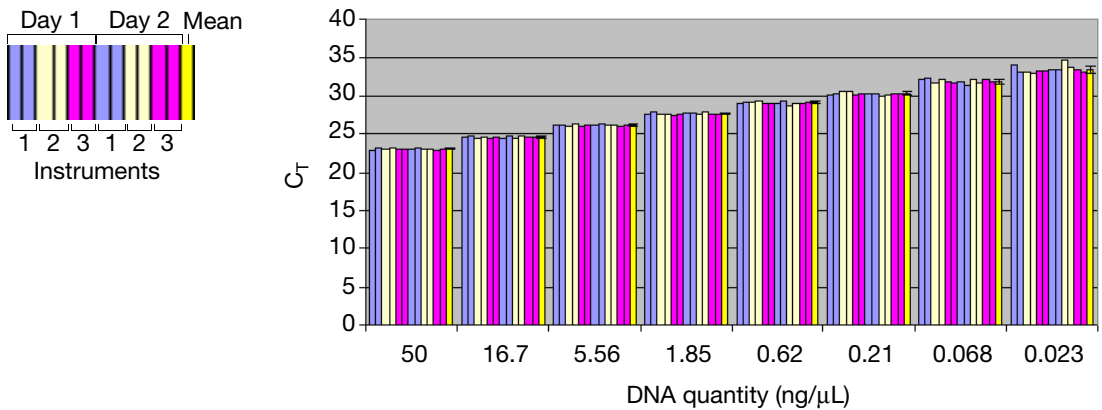


Figure 6-1 Precision using the Quantifiler Human kit

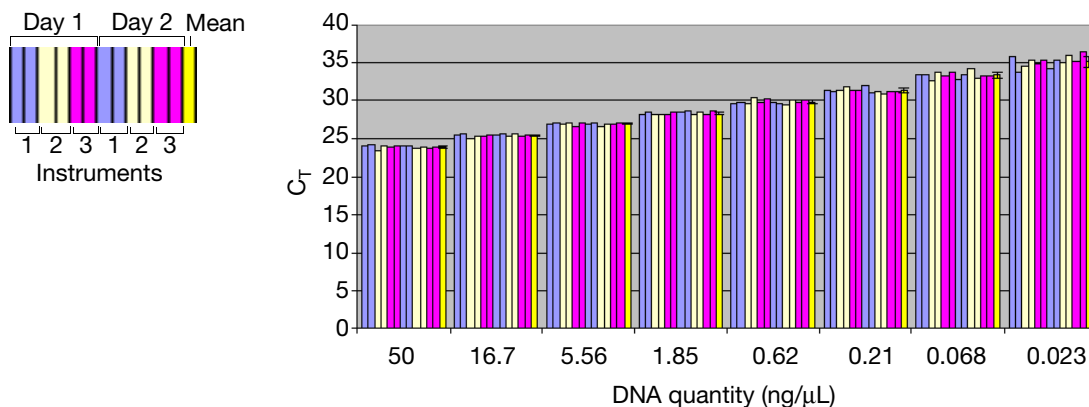


Figure 6-2 Precision using the Quantifiler Y kit

The data show that at lower DNA concentrations, the C_T values increased and the standard deviation increased, most likely because of stochastic effects.

For each sample, the C_T values obtained using the Quantifiler Human kit are lower than those obtained using the Quantifiler Y kit because there are two copies of the autosomal human target locus and only one copy of the Y chromosome target locus.

The C_T values do not vary significantly from run to run or from instrument to instrument. The C_T value from one sample run on three different 7000 instruments varies with an average standard deviation of 0.3. Systematic differences between instruments, which are normally insignificant, are not expected to affect final sample quantification results because, when samples and quantification standards are run on the same plate and instrument, the C_T values are affected equally.

6.1.2 Reproducibility

Experiment Six different human DNA samples were tested for reproducibility of the quantification results.

Table 6-2 Human DNA samples tested for reproducibility

DNA	Sex	Extraction Source
007	Male	Blood
9948	Male	Cell line
Human genomic	Male	Blood
Raji (Lot 1)	Male	Cell line
Raji (Lot 2)	Male	Cell line
K-562	Female	Cell line

Using the concentrations provided by the supplier, the DNA samples were diluted to 2.0 ng/μL (A), 0.5 ng/μL (B), and 0.1 ng/μL (C).

Note: All dilutions were made in T₁₀E_{0.1} Buffer with 20 μg/mL glycogen added as a carrier and stabilizer.

All samples and dilutions were tested in successive runs using the Quantifiler Human kit and the Quantifiler Y kit. Three different runs were performed. Each assay contained two reactions for each of the quantification standards and one reaction for each of the samples.

For each sample reaction the C_{T FAM} values were obtained and the DNA quantity calculated. The mean quantity and standard deviations were calculated for each sample. The 95% confidence interval values were calculated as the mean of the DNA quantity ± two standard deviation units for each sample and expressed as a percentage of the mean quantification result.

Results The following tables show the DNA quantity calculated for all samples and dilutions tested for all three runs using the Quantifiler Human kit (Table 6-3) and the Quantifiler Y kit (Table 6-4).

Table 6-3 Reproducibility using the Quantifiler Human kit

Sample	DNA Quantity (ng/ μ L)				Standard Deviation	95% Confidence (\pm percent)
	Run 1	Run 2	Run 3	Mean		
007 A	2.580	2.830	2.900	2.770	0.168	12.15
007 B	0.894	0.779	0.892	0.855	0.066	15.40
007 C	0.216	0.160	0.192	0.189	0.028	29.68
9948 A	2.300	2.240	2.210	2.250	0.046	4.07
9948 B	0.504	0.481	0.573	0.519	0.048	18.44
9948 C	0.123	0.132	0.132	0.129	0.005	8.06
Human genomic A	1.810	1.790	2.240	1.947	0.254	26.12
Human genomic B	0.495	0.468	0.504	0.489	0.019	7.66
Human genomic C	0.128	0.106	0.106	0.113	0.013	22.41
K-562 A	1.360	1.350	1.360	1.357	0.006	0.85
K-562 B	0.379	0.425	0.460	0.421	0.041	19.28
K-562 C	0.096	0.126	0.096	0.106	0.017	32.42
Raji-1 A	1.920	1.800	1.770	1.830	0.079	8.67
Raji-1 B	0.484	0.402	0.466	0.451	0.043	19.13
Raji-1 C	0.149	0.120	0.104	0.124	0.023	36.69
Raji-2 A	1.720	1.860	1.700	1.760	0.087	9.91
Raji-2 B	0.419	0.407	0.408	0.411	0.007	3.24
Raji-2 C	0.113	0.088	0.061	0.087	0.026	59.50

Table 6-4 Reproducibility using the Quantifiler Y kit

Sample	DNA Quantity (ng/ μ L)				Standard Deviation	95% Confidence
	Run 1	Run 2	Run 3	Mean		
007 A	3.760	3.600	3.840	3.733	0.122	6.55
007 B	1.180	0.898	1.040	1.039	0.141	27.13
007 C	0.238	0.185	0.172	0.198	0.035	35.26
9948 A	2.590	2.540	2.670	2.600	0.066	5.04
9948 B	0.810	0.612	0.709	0.710	0.099	27.88
9948 C	0.146	0.130	0.151	0.142	0.011	15.41
Human genomic A	2.010	1.770	1.760	1.847	0.142	15.33
Human genomic B	0.577	0.462	0.591	0.543	0.071	26.06
Human genomic C	0.081	0.053	0.052	0.062	0.017	54.04
K-562 A	—	—	—	n.d. ^a	n.d.	n.d.
K-562 B	—	—	—	n.d.	n.d.	n.d.
K-562 C	—	—	—	n.d.	n.d.	n.d.
Raji-1 A	2.500	2.090	2.400	2.330	0.214	18.35
Raji-1 B	0.679	0.481	0.565	0.575	0.099	34.57
Raji-1 C	0.123	0.096	0.148	0.122	0.026	42.80
Raji-2 A	2.630	2.050	2.190	2.290	0.303	26.43
Raji-2 B	0.574	0.536	0.612	0.574	0.038	13.24
Raji-2 C	0.091	0.123	0.160	0.125	0.034	55.02

a. n.d. = not determined

The 95% confidence interval shows the approximate range expected for results when using the Quantifiler kits. The average 95% confidence interval for each kit:

- Quantifiler Human kit: $\pm 18.5\%$
- Quantifiler Y kit: $\pm 26.9\%$

The data show that as the DNA concentration decreases, the amount of variability in the quantification results increases. This results from stochastic effects—the statistical principles involved when testing DNA samples with low concentrations. Stochastic effects may cause imbalance or dropouts of alleles when performing STR analysis of DNA samples with low concentrations.

6.1.3 Specificity with a Human DNA Panel

Purified genomic DNA samples from 500 human individuals were obtained from two different commercial sources. Many of the samples were extracted from cell lines that provide distinct genotypes for forensic validation work; other samples were extracted from blood specimens. The sex of all samples was confirmed by genotypic analysis using the AmpF \mathbb{L} STR \mathbb{R} Identifiler \mathbb{R} PCR Amplification Kit (amelogenin locus).

Experiment

Approximately 20 to 40 ng of purified genomic DNA from the Human DNA Panel was used for each Quantifiler kit reaction.

Sequence Detection Systems (SDS) software was used to analyze the data and calculate the $C_{T\text{ FAM}}$ value:

$C_{T\text{ FAM}}$ Value	Result
$C_{T\text{ FAM}} < 40$	+
No amplification after 40 cycles	–

Results The results in Table 6-5 show that:

- The Quantifiler Human kit detected all 500 human DNA samples.
- The Quantifiler Y kit detected all 240 male DNA samples and none of the female DNA samples.

Table 6-5 Specificity with human DNA panel

Sex	Result	
	Quantifiler Human Kit	Quantifiler Y Kit
Male (240)	+	+
Female (260)	+	–

6.1.4 Specificity with a Non-Human Panel

Samples were obtained either as purified DNA or as whole blood from individual animals. For some of the purified DNA samples, the sex of the donor animals was unknown; for remaining samples, the sex and identity of the animals was known. For some species, multiple individuals were tested.

Experiment For many of the reactions, approximately 0.25 to 1.0 ng of DNA was used in each reaction. For a few reactions, up to 40 ng of DNA was used in one reaction.

SDS software was used to analyze the data and calculate the $C_{T\text{ FAM}}$ value:

$C_{T\text{ FAM}}$ Value	Result
$C_{T\text{ FAM}} < 40$	+
No amplification after 40 cycles	–

Results The two human control samples that were tested show expected results (as shown in [Table 6-5 on page 6-11](#)).

Quantifiler Human Kit Results

The Quantifiler Human kit detected DNA from humans and apes, with some less-efficient detection of one other primate. The Quantifiler Human kit:

- Detected DNA from all of the higher ape DNA samples (chimpanzee, gorilla, and orangutan) at an efficiency similar to that of humans
- Detected DNA from macaque monkeys at a significantly reduced efficiency, possibly because of partial homology between the primers and probe and the macaque DNA
- Did not detect DNA from the remaining species

Quantifiler Y Kit Results

The Quantifiler Y kit detected DNA from male humans and chimpanzees but from no other species tested.

Of the DNA samples that were detected using the Quantifiler Human Kit (gorilla, chimpanzee, orangutan, and macaque), the Quantifiler Y kit:

- Detected DNA from the chimpanzees
- Did not detect DNA from the male gorilla
- Did not detect DNA from the female orangutans or macaques

Table 6-6 Specificity with non-human panel

Organism	Sex	Result	
		Quantifiler Human Kit	Quantifiler Y Kit
Gorilla (2)	Female ^a	+	—
Chimpanzee (2)	Unknown	+	+
Orangutan (2)	Female ^a	+	—
Macaque (2)	Female ^a	± ^b	—

Table 6-6 Specificity with non-human panel (*continued*)

Organism	Sex	Result	
		Quantifiler Human Kit	Quantifiler Y Kit
Cat	Unknown	–	–
Dog	Unknown	–	–
Pig	Unknown	–	–
Cow	Unknown	–	–
Mouse	Unknown	–	–
Rabbit	Unknown	–	–
Hamster	Unknown	–	–
Rat	Unknown	–	–
Chicken	Unknown	–	–
Fish	Unknown	–	–
Gorilla	Male	+	–
Cat	Male	–	–
Dog (2)	Male	–	–
Mouse	Male	–	–
Rabbit	Male	–	–
Rat	Male	–	–
Horse (2)	Male	–	–
Bovine	Male	–	–
Sheep	Male	–	–
Pig	Male	–	–
Deer	Male	–	–

Table 6-6 Specificity with non-human panel (*continued*)

Organism	Sex	Result	
		Quantifiler Human Kit	Quantifiler Y Kit
Chicken	Male	–	–
Human	Female	+	–
Human	Male	+	+

a. Sex confirmed by STR analysis.

b. Weak but positive amplification with higher C_T values and lower R_n values than normal for the input amount of DNA in the reaction.

6.1.5 Specificity with a Bacterial Pools Panel

The bacterial pools panel contained purified genomic DNA from 53 bacterial species and one yeast species. The panel included:

- Common gram-negative and gram-positive species
- Species associated with the human gut (for example, *Proteus*, *Providencia*, *Alcaligenes*)
- Species associated with food (*Lactobacillus* spp.)
- Species associated with spoilage and decomposition (for example, *Pseudomonas*, *Flavobacterium*, *Clostridium*, *Candida*)
- Species associated with human enteric disease (for example, *Salmonella*, *Escherichia coli*, *Yersinia*).
- Several species of *Bacillus*, a common and pervasive bacterial genus

Experiment There were approximately 1×10^5 genome copies of each species in each reaction.

SDS software was used to analyze the data and calculate the $C_{T \text{ FAM}}$ value:

$C_{T \text{ FAM}}$ Value	Result
$C_{T \text{ FAM}} < 40$	+
No amplification after 40 cycles	–

Results The Quantifiler Human Kit and the Quantifiler Y kit did not detect DNA from any of the bacterial or yeast species tested.

Table 6-7 Specificity with bacterial pools panel

Species Composition	Result	
	Quantifiler Human Kit	Quantifiler Y Kit
<i>Lactobacillus acidophilus</i> , <i>Lactobacillus delbrueckii</i> (2), <i>Lactobacillus rhamnosus</i> , <i>Lactobacillus casei</i>	—	—
<i>Brochothrix thermosphacta</i> , <i>Brochothrix campestris</i> , <i>Aerococcus viridians</i> , <i>Kurthia gibsonii</i> , <i>Alcaligenes faecalis</i>	—	—
<i>Bacillus subtilis</i> , <i>Bacillus cereus</i> , <i>Bacillus licheniformis</i> , <i>Bacillus mycoides</i> , <i>Bacillus stearothermophilus</i>	—	—
<i>Pseudomonas fluorescens</i> , <i>Flavobacterium odoratum</i> , <i>Clostridium sporogenes</i> , <i>Candida kefyr</i> (yeast), <i>Deinococcus radiodurans</i>	—	—
<i>Lactococcus lactis</i> , <i>Bordetella bronchiseptica</i> , <i>Acinetobacter baumannii</i> , <i>Aeromonas caviae</i> , <i>Corynebacterium varibile</i>	—	—
<i>Nocardia asteroides</i> , <i>Stenotrophomonas maltophilia</i> , <i>Bacillus coagulans</i> , <i>Rhodococcus equi</i> , <i>Acinetobacter calcoaceticus</i>	—	—
<i>Propionibacterium acnes</i> , <i>Clostridium difficile</i> , <i>Fuseobacterium necrophorum</i> , <i>Burkholderia cepacia</i> , <i>Delftia acidovorans</i>	—	—
<i>Micrococcus luteus</i> , <i>Streptomyces rimosus</i> , <i>Gordonia sputi</i> , <i>Legionella ansia</i> , <i>Pasteurella aerogenes</i>	—	—
<i>Citrobacter freundii</i> , <i>Klebsiella pneumoniae</i> , <i>Escherichia hermanii</i> , <i>Enterobacter cloacae</i> , <i>Escherichia coli</i> O157:H7	—	—
<i>Salmonella enteritidis</i> , <i>Shigella dysenteriae</i> , <i>Proteus vulgaris</i> , <i>Pseudomonas aeruginosa</i> , <i>Hafnia alvei</i>	—	—
<i>Yersinia enterocolitica</i> , <i>Campylobacter coli</i> , <i>Providencia stuartii</i> , <i>Vibrio parahaemolyticus</i> , <i>Alcaligenes faecalis</i>	—	—

6.1.6 Sensitivity

Human genomic DNA samples were obtained from different commercial sources. For each DNA sample, a dilution series was made and each dilution was tested with the Quantifiler Human kit and the Quantifiler Y kit.

DNA Samples Tested

Five different human DNA samples were tested.

Table 6-8 Human DNA samples tested for sensitivity

Sample	Extraction Source
007	Human male blood
9948	Human male cell line
Human genomic	Human male blood
Raji	Human male cell line
K-562	Human female cell line

Experiment

Using the concentrations provided by the suppliers, five-fold serial dilutions of the DNA samples were made. Concentrations ranged from 10 ng/ μ L to 0.016 ng/ μ L (16 pg/ μ L).

Note: All dilutions were made in $T_{10}E_{0.1}$ Buffer with 20 μ g/mL glycogen added as a carrier and stabilizer.

For each 25- μ L reaction, 2.0 μ L of DNA sample was used.

Results

A plot of the C_T values versus the known DNA quantities showed the expected log-linear relationship between the two quantities. All dilutions, including samples at the lowest concentration (16 pg/ μ L), gave positive results for the Quantifiler Human kit and the Quantifiler Y kit. For each dilution series, the data points formed an acceptable standard curve. The small differences in C_T values among the dilutions of different DNA samples likely reflect differences in the quantification measurements made by each supplier.

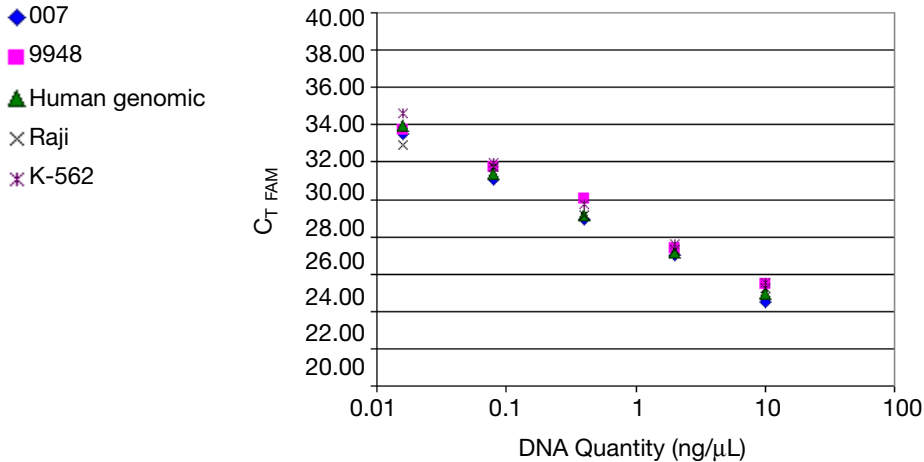


Figure 6-3 Sensitivity using the Quantifiler Human kit

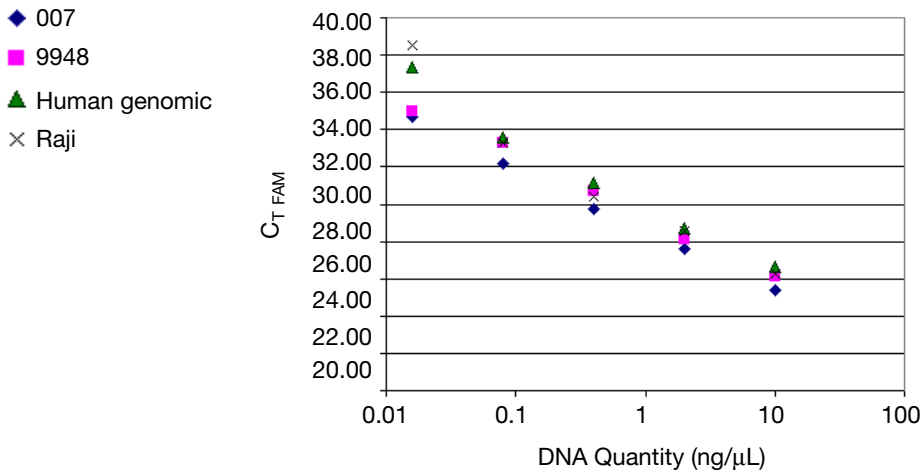


Figure 6-4 Sensitivity using the Quantifiler Y kit

6.1.7 Stability

DNA samples from various origins are commonly contaminated with organic and inorganic compounds that inhibit the amplification of nucleic acids by PCR. These PCR inhibitors can interfere with the reaction and cause varying levels of reduced PCR efficiency, including complete inhibition of PCR. A wide variety of PCR inhibitors has been reported, including in DNA samples extracted

from blood stains. One example is hematin, which has been found in DNA samples extracted from blood stains. Because the solubility of hematin is similar to that of DNA, it is thought that it is extracted and purified with the DNA. The presence of hematin in DNA samples may interfere with PCR by inhibiting polymerase activity.

Bovine serum albumin (BSA) is used in enzymatic reactions because it appears to increase the efficiency of the PCR reaction, most likely acting as a chelating agent with many inhibitors. BSA is added to the Quantifiler kit and AmpF Φ STR[®] kit reaction mixes specifically to counteract the presence of PCR inhibitors.

Experiment Human genomic DNA was mixed with varying concentrations of hematin: 0 μ M, 10 μ M, 12 μ M, 14 μ M, 16 μ M, 18 μ M, 20 μ M, and 40 μ M. 2.0 μ L of each DNA/hematin mix, containing 1.0 ng total of human DNA, was quantified using the Quantifiler Human kit and Quantifiler Y kit; the same amounts of samples were added to reactions using the AmpF Φ STR[®] Identifiler[®] PCR Amplification Kit. Identifiler kit reactions were analyzed on a 3100 instrument. Data were analyzed with GeneScan[®] Software v3.7.1 and Genotyper[®] Software v3.7, for use with the Windows NT[®] operating system.

Results Amplification plots (Figures 6-5 and 6-6) showed lower ΔR_n values and higher C_T values as the concentration of hematin increased. C_T results and corresponding quantification results were relatively stable up to 14 μ M hematin, with results more affected at higher concentrations. As the concentration of hematin increased, the PCR efficiency in the Quantifiler kit reactions and the Identifiler kit reactions decreased. For the Quantifiler Human kit, complete inhibition occurred at 40 μ M, and for the Quantifiler Y kit, complete inhibition occurred at 18 μ M, 20 μ M, and 40 μ M. The inhibition may be stronger with the Quantifiler Y kit because there is only one copy of the haploid Y chromosome target locus for the Quantifiler Y kit and two copies of the diploid autosomal target locus for the Quantifiler Human kit.

The IPC system is more sensitive to PCR inhibition. For the Quantifiler Human kit, in samples containing more than 16 μ M hematin, amplification of IPC detectors failed. In samples containing less hematin, amplification of IPC detectors was inhibited (Figure 6-7). Although the Human detector amplified for the 16 μ M, 18 μ M and 20 μ M hematin samples, the failure of IPC amplification

in those reactions indicates that the presence of PCR inhibitors is likely. Because the IPC system components are the same in both Quantifiler kits, the IPC results for the Quantifiler Y kit were similar to those for the Quantifiler Human kit.

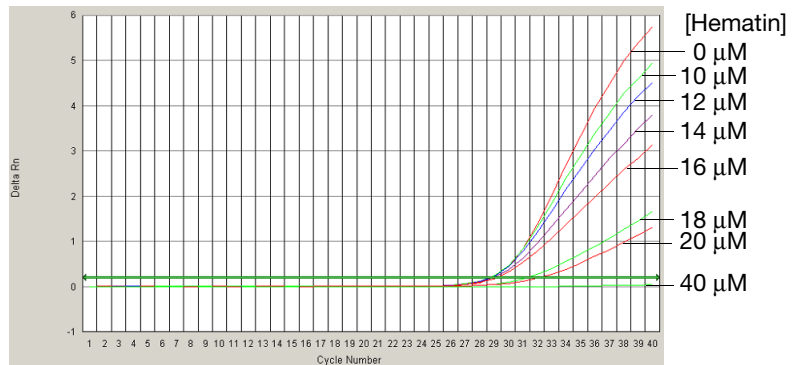


Figure 6-5 Inhibition studies: Quantifiler Human kit

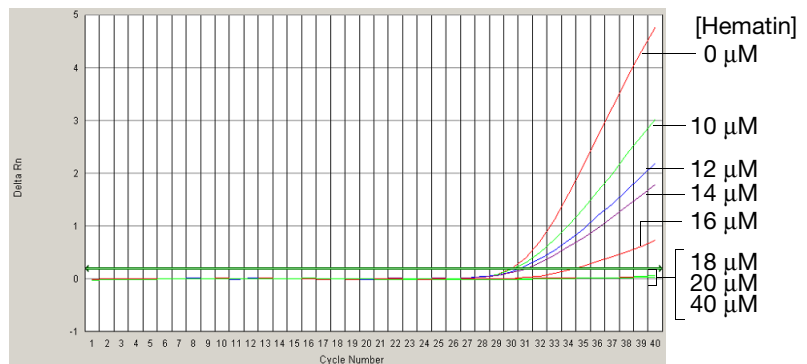


Figure 6-6 Inhibition studies: Quantifiler Y kit

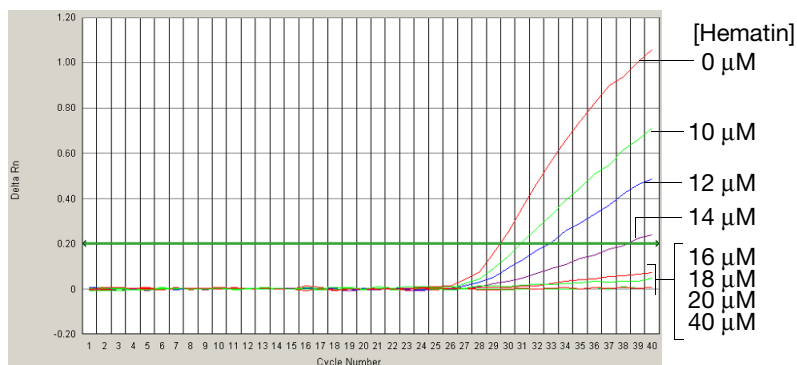


Figure 6-7 Inhibition studies: IPC detector

The results of STR analysis using the Identifiler kit ([Figure 6-8](#)) were consistent with the results from the Quantifiler kits: as the concentration of hematin increased, the overall STR peak profile decreased. Complete STR profiles were obtained at hematin concentrations up to 20 μM . The STR amplification reaction was completely inhibited by 40 μM hematin. The results from the Quantifiler kits provided reasonable predictions of samples that would fail STR analysis because of the presence of the PCR inhibitor. The STR profiles for the positive and negative controls are included for reference.

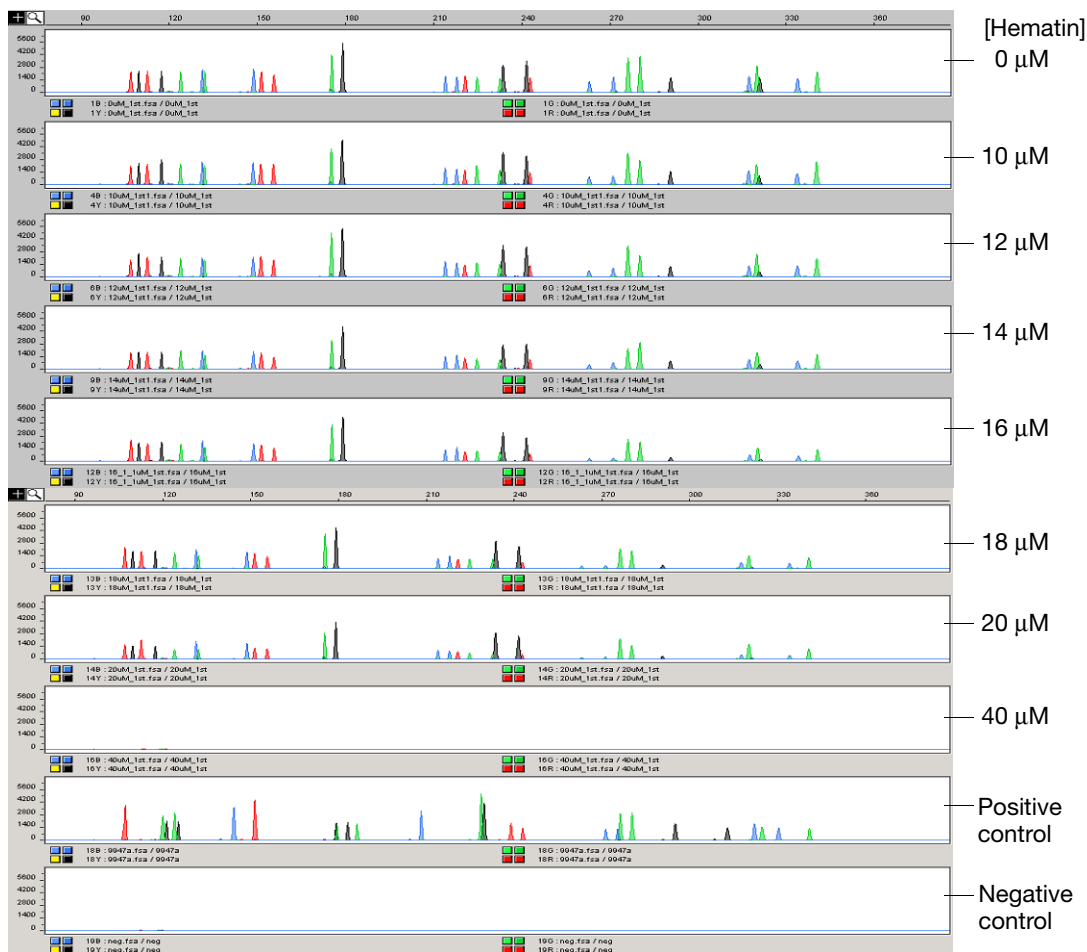


Figure 6-8 Inhibition studies: STR analysis

6.1.8 Mixture Studies

The mixture studies in this section were designed to simulate circumstances in which a small component of male DNA must be discerned from a high background of female DNA. Evidence samples may contain DNA from more than one individual, and this should be considered when interpreting the results. Applied Biosystems recommends that individual laboratories assign a minimum peak height threshold based on validation experiments performed in each laboratory.

Experiment Purified genomic DNA from the Raji (male) and K-562 (female) cell lines were mixed in ratios of 1:1, 1:4, 1:16, 1:64, 1:256 and 1:1024 (Raji:K-562). The male DNA was added at a constant level of 0.05 ng/μL in all samples, and the female DNA was present at amounts ranging from 0.05 ng/μL in the 1:1 sample to 50 ng/μL in the 1:1024 sample. The DNA amounts were calculated based only on the DNA concentrations provided by the suppliers and were not calibrated with the Quantifiler kits.

The mixtures were tested with the Quantifiler Human kit and the Quantifiler Y kit to determine the concentrations of total human genomic DNA (Quantifiler Human kit) and male DNA only (Quantifiler Y kit). For each sample, three replicate reactions were performed for each assay. Each assay used the same set of 8 human genomic DNA quantification standards run in duplicate reactions for each assay and both assays were run on the same reaction plate. The reaction plates were run on a 7000 instrument.

Results The quantification results ([Figure 6-9 on page 6-23](#)) from using the Quantifiler Human kit varied from an average of 0.16 ng/μL for the 1:1 sample to 38 ng/μL for the 1:1024 sample, consistent with the increasing amounts of female DNA present.

The quantification results from using the Quantifiler Y kit varied from between 0.034 ng/μL to 0.063 ng/μL for all samples, regardless of the amount of female DNA present.

For the 1:1024 sample, the results showed a ratio of male DNA to total DNA of 1:974. Differences between target concentrations and actual measurements were expected because the amounts of DNA added to the mixtures were based only on the DNA concentrations provided by the suppliers and were not calibrated with the Quantifiler kits.

In all samples, the male DNA was detected and quantified accurately, regardless of the amount of female DNA present.

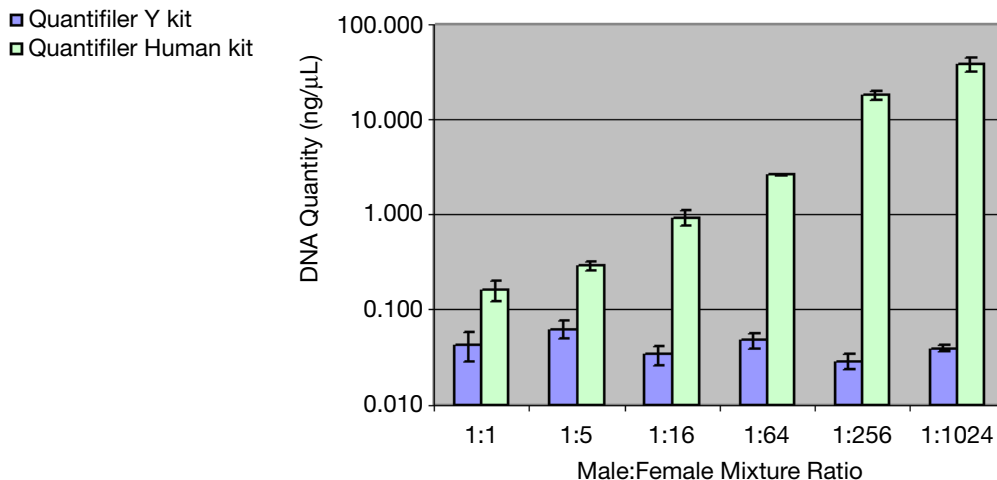


Figure 6-9 DNA quantities determined in mixture studies

6.1.9 Degraded DNA Studies

Forensic samples may be exposed to environmental conditions that degrade DNA molecules and reduce their amplification efficiency in PCR reactions. Exposure to environmental conditions can reduce the overall DNA concentration and may cause fragmentation of full-length DNA molecules into smaller fragments. DNA fragmentation makes it difficult to amplify longer segments such as the larger STR loci. Because of such potential occurrences, the validation of forensic DNA methods often involves studies of the effects of degradation on the amplification and detection of DNA.

The Quantifiler® kits were tested with DNA degraded with the DNA nuclease DNase I. The degraded DNA samples were tested with the Quantifiler Human kit and the Quantifiler Y kit to determine the quantity of amplifiable DNA in each time point. Results obtained using the Quantifiler kits were used to calculate DNA input for an STR assay using an ABI PRISM® 3100 Genetic Analyzer.

Experiment

A time-course of exposure to DNase I was performed on a sample of high molecular weight human genomic DNA to generate a series of samples with varying levels of degradation. The time points in the DNase I treatment were 0 minutes (untreated), 1 minute, 2 minutes, 3 minutes, 4 minutes, 5 minutes, 15 minutes and 60 minutes. Samples from all time points were run on a 2% agarose gel for 25 minutes and visualized by staining with ethidium bromide. The

treated DNA samples were examined by agarose gel electrophoresis to determine the average size of the DNA fragments at each time point. The degraded DNA samples were tested with the Quantifiler Human kit and the Quantifiler Y kit to determine the quantity of amplifiable DNA in each time point.

Using the results from the Quantifiler kits, the volumes of DNA required for AmpF ℓ STR $^{\text{®}}$ Identifier $^{\text{®}}$ kit reactions were calculated so that 1.0 ng/ μ L was added for each reaction. The PCR products were run on an ABI PRISM $^{\text{®}}$ 3100 Genetic Analyzer.

Results Agarose gel electrophoresis showed that the DNase I treatment reduced the average size of DNA fragments to 100 basepairs (bp) or less within the first 5 minutes (Figure 6-10).

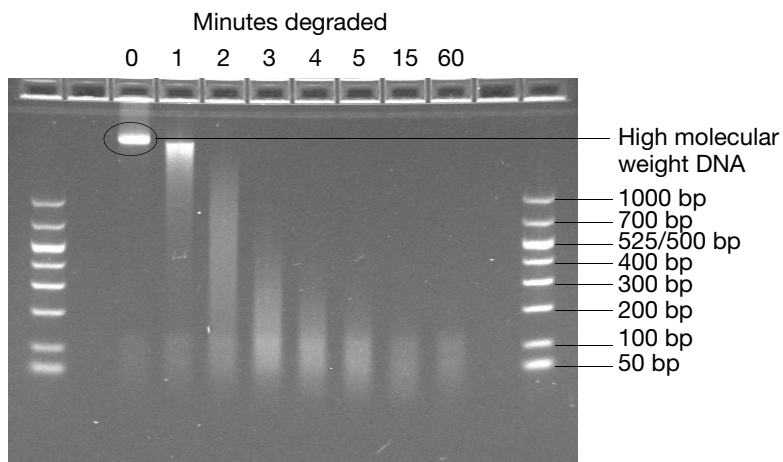


Figure 6-10 DNase I degradation of human genomic DNA

The results from the Quantifiler kits (Figures 6-11 and 6-12) showed higher C_T values with longer DNase exposure times, corresponding to lower amounts of amplifiable DNA in the samples. According to results from the Quantifiler Human kit, the amount of amplifiable DNA decreased from 12.0 ng/ μ L to 1.2 ng/ μ L at the 5-minute time point and to 0.11 ng/ μ L at the 15-minute time point. At the 60-minute time point, no amplifiable DNA was detected.

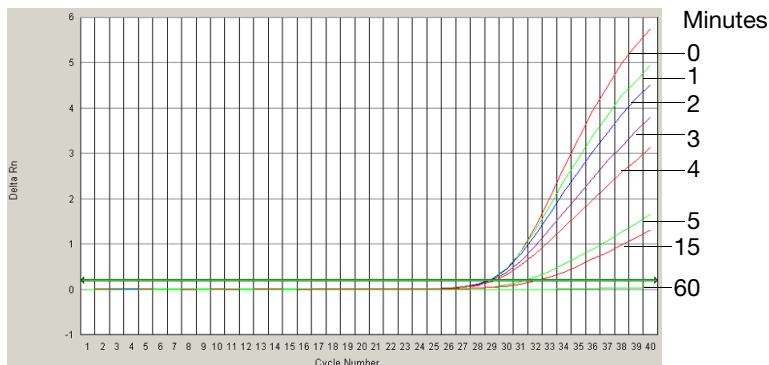


Figure 6-11 Degraded DNA: Quantifiler Human kit amplification plot

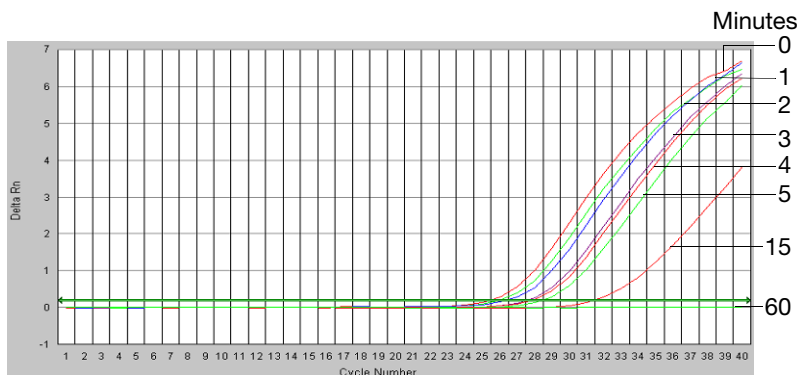


Figure 6-12 Degraded DNA: Quantifiler Y kit amplification plot

Using the DNA quantification results from the Quantifiler Human kit, 1.0 ng of each DNA sample was added to Identifiler kit reactions. As the concentration of amplifiable DNA decreased because of degradation, the sample volume required in the reaction increased.

Identifiler kit results at 1.0 ng/μL produced complete STR profiles up to the 5-minute time point, although the amount of amplifiable DNA (according to the Quantifiler kit) was reduced by 90% relative to the untreated control ([Figure 6-13](#)). The peak heights were reduced for the more degraded samples, but profiles were still detected. The 15-minute time point contained only 1% of the original amount of

amplifiable DNA and produced only a partial STR profile of mostly smaller molecular weight loci. At 60 minutes, no DNA was detected by the Quantifiler kits (Figures 6-11 and 6-12) or the Identifiler kit (Figure 6-13).

The Quantifiler kits can be used to report the amount of amplifiable DNA in a sample but not the amount of DNA degradation. Using the quantification data from the kits to determine the amount of sample input for STR analysis may help to correct for the loss of amplifiable DNA because of degradation, but if the level of DNA degradation is so high that the remaining DNA fragments are too small, the sample will not amplify by using the Quantifiler kits or the STR kits.

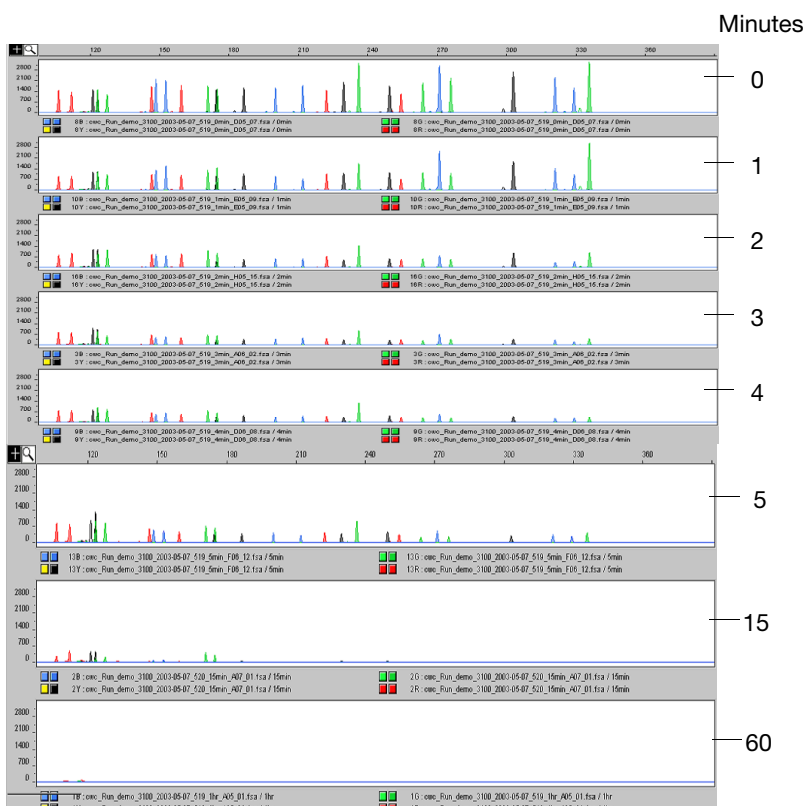


Figure 6-13 STR analysis using degraded DNA

6.1.10 Comparisons with Other Methods

Purified DNA samples were quantified using the Quantifiler Human kit and the Quantifiler Y kit. The results were compared to results obtained from measuring absorbance at 260 nm (A_{260}), using a dye intercalation method, and using the Quantiblot® Human DNA Quantitation Kit (Applied Biosystems).

The methods tested show different sensitivity ranges and different specificities.

Table 6-9 Comparison: sensitivity and specificity of methods

Method	Sensitivity	Specificity
A_{260}	Cannot detect DNA in the picogram range.	Not specific for human genomic DNA. Detects single-stranded DNA, double-stranded DNA, and RNA.
Dye intercalation	25 pg/mL ^a	Not specific for human genomic DNA
Quantiblot kit	2 ng/μL to 0.03125 ng/μL	Specific for human genomic DNA

a. Value obtained from the manufacturer's documentation.

6.1.11 Comparison with A_{260} and Quantiblot Kit

The concentration of DNA was measured for 50 human genomic DNA samples using a A_{260} method, the Quantiblot kit, and the Quantifiler kits. The DNA quantification results were compared.

Resolution Panel The resolution panel, a set of 50 human genomic DNA samples purified from blood, was tested. The samples were database type samples because they were extracted from blood specimens and had uniform high concentrations of DNA between approximately 10 and 20 ng/μL. All samples were within the range of sensitivity for the A_{260} method.

Experiment Each DNA sample was quantified using:

- **A₂₆₀ method** – Absorbance at 260 nm was measured. DNA concentration was calculated using the formula:

$$\text{Concentration } (\mu\text{g/mL}) = 50 \times A_{260}$$
- **Quantiblot kit** – DNA was quantified using a protocol for chemiluminescence detection with film autoradiography.
- **Quantifiler kits** – DNA was quantified using the standard procedure.

For each sample, the percent differences between Quantifiler kit results and results from the other two methods were calculated. The differences were expressed as a percentage of the reference method. For each method, the average percent differences from Quantifiler kit results were calculated. For comparisons with the Quantifiler Y kit, only results from male samples were used.

Results Table 6-10 shows the DNA quantification results for all 50 samples in the resolution panel and for the three methods. The table also shows the percent differences between the results from the Quantifiler kits and the other two methods. There is no A₂₆₀ data for two samples (13 and 17), and all female samples were excluded from the comparisons to the Quantifiler Y kit results.

Table 6-10 Comparison with A₂₆₀ and Quantiblot kit

Sample	Sex	A ₂₆₀ Result (ng/μL)	QB ^a Result (ng/μL)	Quantifiler Human Kit			Quantifiler Y Kit		
				Result (ng/μL)	% Diff. from A ₂₆₀	% Diff. from QB	Result (ng/μL)	% Diff. from A ₂₆₀	% Diff. from QB
1	M	17.5	20	6.69	61.7	66.6	10.13	41.9	49.4
2	M	15.4	20	14.3	7.1	28.5	16.78	9.0	16.1
3	M	13.9	30	15.48	11.4	48.4	14.30	2.9	52.3
4	M	11.4	20	12.44	9.6	37.8	12.45	9.7	37.8
5	M	10.3	20	12.69	23.2	36.6	11.00	6.8	45.0
6	M	13.9	20	12.54	9.8	37.3	13.56	2.4	32.2
7	M	11.5	40	13.78	20.1	65.6	12.28	7.1	69.3

Table 6-10 Comparison with A₂₆₀ and Quantiblot kit (*continued*)

Sample	Sex	A ₂₆₀ Result (ng/μL)	QB ^a Result (ng/μL)	Quantifiler Human Kit			Quantifiler Y Kit		
				Result (ng/μL)	% Diff. from A ₂₆₀	% Diff. from QB	Result (ng/μL)	% Diff. from A ₂₆₀	% Diff. from QB
8	M	11.2	20	13.51	21.2	32.5	11.77	5.6	41.2
9	M	9.8	20	15.09	54.0	24.6	13.06	33.3	34.7
10	M	9.7	20	13.98	44.1	30.1	12.29	26.7	38.6
11	M	13.0	20	11.27	13.3	43.7	12.85	1.2	35.8
12	M	13.3	30	9.92	25.1	66.9	11.59	12.5	61.4
13	M	nd	14	13.90	n.d.	0.7	11.31	n.d.	19.2
14	F	12.8	16	13.90	9.0	13.1	neg	n.d.	n.d.
15	M	15.7	16	12.62	19.4	21.1	13.89	11.2	13.2
16	M	12.1	24	13.09	8.2	45.5	10.78	10.9	55.1
17	M	nd	20	12.81	n.d.	36.0	14.36	n.d.	28.2
18	M	13.5	24	8.18	39.4	65.9	10.25	24.1	57.3
19	M	13.2	20	10.37	21.4	48.2	13.12	0.6	34.4
20	M	12.9	16	12.69	1.2	20.7	12.36	3.8	22.8
21	M	11.0	14	13.48	22.9	3.7	13.00	18.5	7.1
22	M	11.5	24	12.23	6.6	49.0	12.85	12.0	46.5
23	M	10.9	14	10.91	0.6	22.1	11.73	8.1	16.2
24	M	12.4	20	15.19	22.8	24.1	14.38	16.2	28.1
25	M	10.8	20	15.21	41.5	24.0	18.07	68.1	9.7
26	F	13.9	20	14.00	1.1	30.0	— ^b	n.d.	n.d.
27	F	11.5	32	13.16	14.4	58.9	—	n.d.	n.d.
28	F	11.5	40	10.51	8.6	73.7	—	n.d.	n.d.

Table 6-10 Comparison with A₂₆₀ and Quantiblot kit (*continued*)

Sample	Sex	A ₂₆₀ Result (ng/μL)	QB ^a Result (ng/μL)	Quantifiler Human Kit			Quantifiler Y Kit		
				Result (ng/μL)	% Diff. from A ₂₆₀	% Diff. from QB	Result (ng/μL)	% Diff. from A ₂₆₀	% Diff. from QB
29	F	11.2	20	10.45	6.3	47.8	–	n.d.	n.d.
30	F	16.0	20	12.56	21.5	37.2	–	n.d.	n.d.
31	F	10.9	20	12.12	11.7	39.4	–	n.d.	n.d.
32	F	10.9	40	9.42	13.6	76.5	–	n.d.	n.d.
33	F	11.5	20	13.95	21.3	30.3	–	n.d.	n.d.
34	F	10.4	20	12.14	16.7	39.3	–	n.d.	n.d.
35	F	11.1	40	12.38	11.3	69.1	–	n.d.	n.d.
36	F	10.5	20	13.38	28.0	33.1	–	n.d.	n.d.
37	F	12.0	24	12.50	4.2	47.9	–	n.d.	n.d.
38	F	10.8	20	9.59	11.0	52.1	–	n.d.	n.d.
39	F	11.4	16	10.42	8.8	34.9	–	n.d.	n.d.
40	F	10.4	40	11.16	7.3	72.1	–	n.d.	n.d.
41	F	12.6	20	12.49	0.9	37.6	–	n.d.	n.d.
42	F	12.5	28	8.68	30.3	69.0	–	n.d.	n.d.
43	F	12.2	20	13.57	11.5	32.2	–	n.d.	n.d.
44	F	9.8	16	9.42	3.9	41.1	–	n.d.	n.d.
45	F	12.4	16	10.96	11.6	31.5	–	n.d.	n.d.
46	F	12.2	16	11.49	5.4	28.2	–	n.d.	n.d.
47	F	10.4	40	12.93	24.1	67.7	–	n.d.	n.d.
48	F	12.3	20	12.23	0.6	38.9	–	n.d.	n.d.

Table 6-10 Comparison with A₂₆₀ and Quantiblot kit (*continued*)

Sample	Sex	A ₂₆₀ Result (ng/μL)	QB ^a Result (ng/μL)	Quantifiler Human Kit			Quantifiler Y Kit		
				Result (ng/μL)	% Diff. from A ₂₆₀	% Diff. from QB	Result (ng/μL)	% Diff. from A ₂₆₀	% Diff. from QB
49	F	10.7	40	15.02	40.4	62.5	–	n.d.	n.d.
50	F	12.8	32	13.50	5.5	57.8	–	n.d.	n.d.

a. Quantiblot kit method

b. Negative (–) result

The different methods produced similar quantification results.

Table 6-11 Average differences from A₂₆₀ and Quantiblot kit

Method	Average Difference (%)	
	Quantifiler Human Kit	Quantifiler Y Kit
A ₂₆₀	16.9	15.1
Quantiblot	42.0	35.5

6.1.12 Comparison with A₂₆₀ and Dye Intercalation

The concentration of DNA was measured for 13 human genomic DNA samples using the A₂₆₀ method, a dye intercalation method, and the Quantifiler kits.

DNA Samples Tested

Six human genomic DNA samples were obtained from different commercial sources.

Table 6-12 Human DNA samples tested with A₂₆₀ and dye intercalation

DNA	Sex	Extraction Source
007	Male	Blood
9948	Male	Cell line

Table 6-12 Human DNA samples tested with A_{260} and dye intercalation

DNA	Sex	Extraction Source
Human genomic	Male	Blood
Raji-1	Male	Cell line
Raji-2	Male	Cell line
K-562	Female	Cell line

Experiment

Using the concentrations provided by the supplier, the DNA samples were diluted to 2.0 ng/ μ L (A), 0.5 ng/ μ L (B), and 0.1 ng/ μ L (C).

Note: All dilutions were made in $T_{10}E_{0.1}$ Buffer with 20 μ g/mL glycogen added as a carrier and stabilizer.

All sample dilutions were quantified using the following methods:

- **A_{260}** – Because the concentrations of the dilutions extended below the detection limit of the spectrophotometer, ultraviolet absorbance at 260 nm was measured for only the highest dilution (2.0 ng/ μ L).

DNA concentration was calculated from the formula:

$$\text{Concentration } (\mu\text{g/mL}) = 50 \times A_{260}$$

The results calculated for the 2.0 ng/ μ L dilutions were then extrapolated for the other dilutions (0.5 ng/ μ L and 0.1 ng/ μ L), using the known dilution factors.

- **Dye intercalation** – The microplate assay mode was used and the plate was read on an ABI PRISM[®] 7700 Sequence Detection System (7700 SDS). All of the sample dilutions were within the detection range of the assay. The assay was run using the λ bacteriophage DNA quantification standard supplied with the kit and a quantification standard based on Raji human genomic DNA. There were significant differences between the standard curves from the λ DNA and Raji DNA. The results obtained from using the Raji DNA standard were used in this experiment because the Raji DNA was considered to be more similar to the DNA measured in these experiments and because the results from using the Raji DNA standard were closer to the results obtained by the other methods.

- **Quantifiler kits** – DNA was quantified using the standard procedure. The Quantifiler® Human DNA standard provided with the kits was used as recommended, with duplicate reactions for each of eight serial dilutions.

For each sample, the percent differences between Quantifiler kit results and results from the other two methods were calculated. The differences were expressed as a percentage of the reference method. For each method, the average percent differences from Quantifiler kit results were calculated. For comparisons with the Quantifiler Y kit, only results from male samples were used.

Results Table 6-13 shows the DNA concentrations calculated for all samples using the A_{260} method, the dye intercalation method, Quantifiler Human kit, and Quantifiler Y kit. It also shows the percent differences calculated for the comparisons between the Quantifiler Human kit or the Quantifiler Y kit and the A_{260} method and the dye intercalation method.

Table 6-13 Comparison with A₂₆₀ and dye intercalation

Sample	A ₂₆₀ Result (ng/μL)	DI ^a Result (ng/μL)	Quantifiler Human Kit			Quantifiler Y Kit		
			Result (ng/μL)	% Diff. from A ₂₆₀	% Diff. from DI	Result (ng/μL)	% Diff. from A ₂₆₀	% Diff. from DI
007 A	2.74	2.502	2.580	5.8	3.1	3.760	37.2	50.3
007 B	0.685	0.756	0.894	30.5	18.3	1.180	72.3	56.2
007 C	0.137	0.176	0.216	57.7	22.6	0.238	73.7	35.1
9948 A	1.9	2.286	2.300	21.1	0.6	2.590	36.3	13.3
9948 B	0.475	0.496	0.504	6.1	1.5	0.810	70.5	63.2
9948 C	0.095	0.103	0.123	29.5	19.4	0.146	53.7	41.7
Human genomic A	2.2	2.270	1.810	17.7	20.3	2.010	8.6	11.5
Human genomic B	0.55	0.584	0.495	10.0	15.2	0.577	4.9	1.1
Human genomic C	0.11	0.134	0.128	16.4	4.8	0.081	26.2	39.6
Raji-1 A	2	1.271	1.920	4.0	51.0	2.500	25.0	96.7
Raji-1 B	0.5	0.351	0.484	3.2	38.1	0.679	35.8	93.7
Raji-1 C	0.1	0.085	0.149	49.0	76.1	0.123	23.0	45.4
Raji-2 A	1.98	1.262	1.720	13.1	36.3	2.630	32.8	108.4
Raji-2 B	0.495	0.357	0.419	15.4	17.3	0.574	16.0	60.7
Raji-2 C	0.099	0.110	0.113	14.1	2.5	0.091	7.7	17.1
K-562 A	2.76	1.317	1.360	50.7	3.3	neg	n.d.	n.d.
K-562 B	0.69	0.365	0.379	45.1	3.9	neg	n.d.	n.d.
K-562 C	0.138	0.104	0.096	30.4	7.9	neg	n.d.	n.d.

a. Dye intercalation method

The different methods produced similar quantification results.

Table 6-14 Average differences from A_{260} and dye intercalation

Method	Average Difference (%)	
	Quantifiler Human Kit	Quantifiler Y Kit
A_{260}	23.3	34.9
Dye intercalation	19.0	48.0

6.1.13 Assay Background

An experiment was performed to check the assay system for false-positive results that would indicate the presence of human DNA in a sample that contained none.

Experiment For each Quantifiler kit, 48 negative control reactions were set up. PCR Mixes were prepared and dispensed into wells of the reaction plate according to the standard procedure. For each negative control reaction, 2 μL of $T_{10}E_{0.1}$ Buffer was added. All standard assay parameters were used, except that the number of thermal cycles was extended from 40 to 50 for increased stringency.

Results [Figures 6-14](#) and [6-15](#) show that all 48 reactions with each assay were negative for their respective human DNA targets. The IPC reactions amplified for all reactions in both assays, indicating that the assay systems performed normally. These data show that there is no inherent false-positive background associated with the Quantifiler kits. However, the assays are extremely sensitive, and achieving clean results requires care in assay setup and good contamination control for reagents, instruments, and laboratory work surfaces.

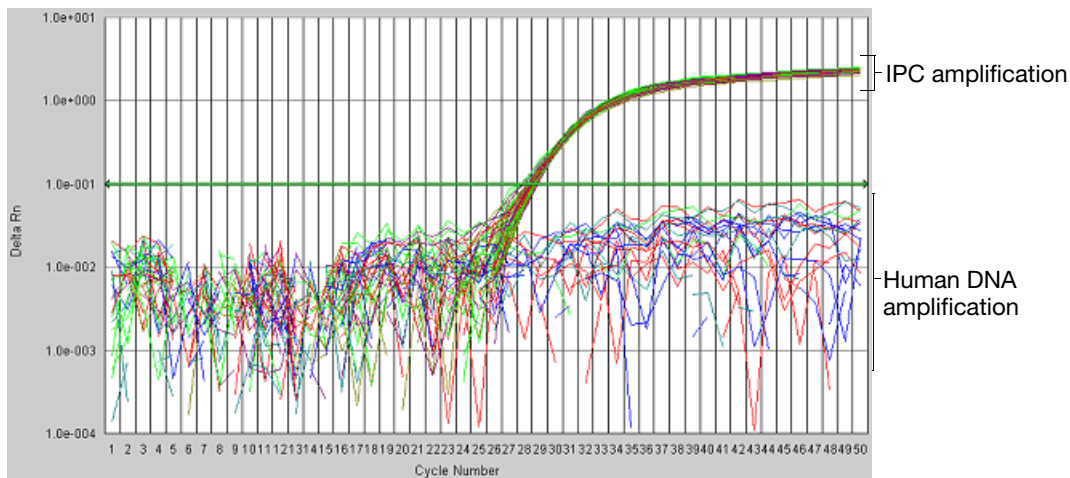


Figure 6-14 Assay background with the Quantifiler Human kit

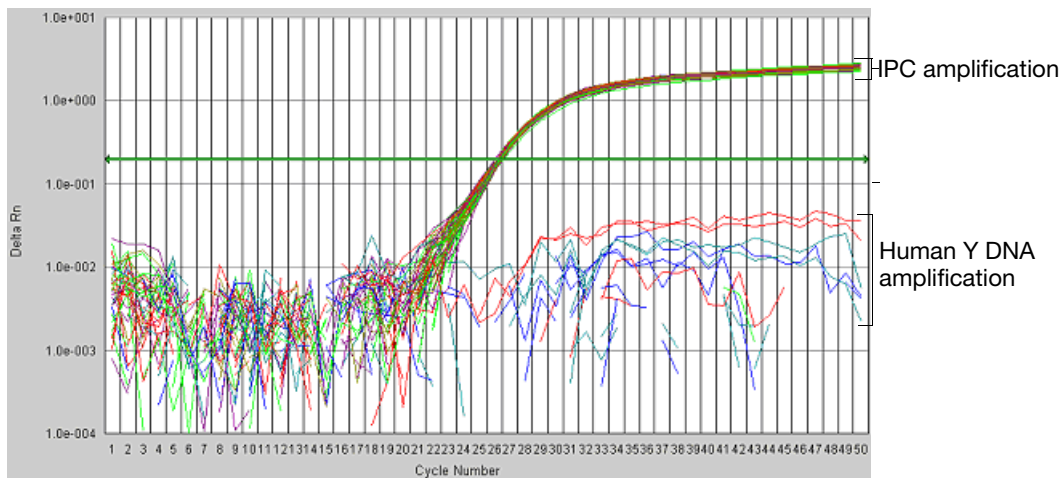


Figure 6-15 Assay background with the Quantifiler Y kit

Section 6.2 Applied Biosystems

7900HT Real-Time PCR System Validation (SDS Software v2.0)

Overview Certain performance parameters for the Quantifiler kits were also tested separately using the ABI PRISM® 7900HT Sequence Detection System (7900HT SDS). The experiments performed for the 7900HT SDS were less exhaustive than those for the 7000 instrument (see previous section) and were performed to test and compare the most sensitive parameters of assay performance between the two instrument platforms.

Experiments	6.2.1 Precision (7900HT SDS)	6-37
	6.2.2 Mixture Studies (7900HT SDS)	6-40
	6.2.3 Comparisons with Other Methods (7900HT SDS)	6-42

6.2.1 Precision (7900HT SDS)

Experiment One set of eight serial dilutions of the Quantifiler Human DNA Standard was prepared. The dilutions ranged from 50 ng/μL to 23 pg/μL in three-fold increments.

Three identical runs containing both Quantifiler Human and Quantifiler Y kits were performed, each containing duplicate reactions of the dilutions for each assay. The three runs were performed on different days on the same 7900HT SDS instrument, all using standard thermal cycler conditions for the Quantifiler kits.

The $C_{T\text{ FAM}}$ values were recorded and the means and standard deviations of the $C_{T\text{ FAM}}$ values were calculated for each of the eight dilutions using the Quantifiler Human kit and the Quantifiler Y kit.

Results Table 6-15 shows the means and standard deviations of the $C_{T\text{ FAM}}$ values calculated for all reactions of each quantification standard dilution for the Quantifiler Human kit and the Quantifiler Y kit.

Table 6-15 Means and standard deviations for C_T results

DNA Quantity (ng/ μ L)	Quantifiler Human Kit		Quantifiler Y Kit	
	C_T (Mean)	Standard Deviation	C_T (Mean)	Standard Deviation
50	23.83	0.13	24.50	0.09
16.7	25.36	0.08	26.08	0.09
5.56	26.79	0.08	27.50	0.06
1.85	28.14	0.08	29.03	0.08
0.62	29.56	0.14	30.68	0.30
0.21	31.00	0.06	32.54	0.42
0.068	32.51	0.25	34.41	0.56
0.023	33.86	0.49	35.59	0.58

The following results are consistent with the 7000 SDS results:

- C_T vs. sample concentration
- Standard deviations of the C_T values
- C_T value calculated using the Quantifiler Human kit was lower than that for the Quantifiler Y kit because there is only one copy of the Y chromosome target locus and two copies of the autosomal human target locus.

The C_T results for all quantification standard dilutions reactions using the Quantifiler Human kit and the Quantifiler Y kit are displayed in Figures 6-16 and 6-17. For each of the dilutions, the mean and the standard deviation of $C_{T\text{FAM}}$ for the repeated runs is shown.

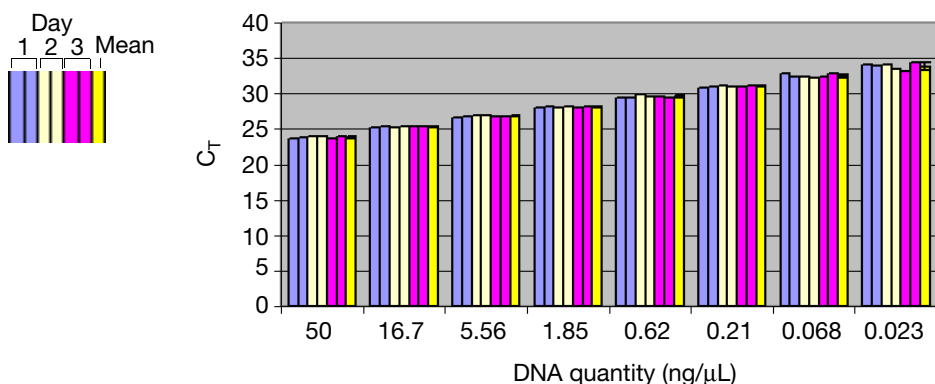


Figure 6-16 Precision: Quantifiler Human kit C_T results (7900HT SDS)

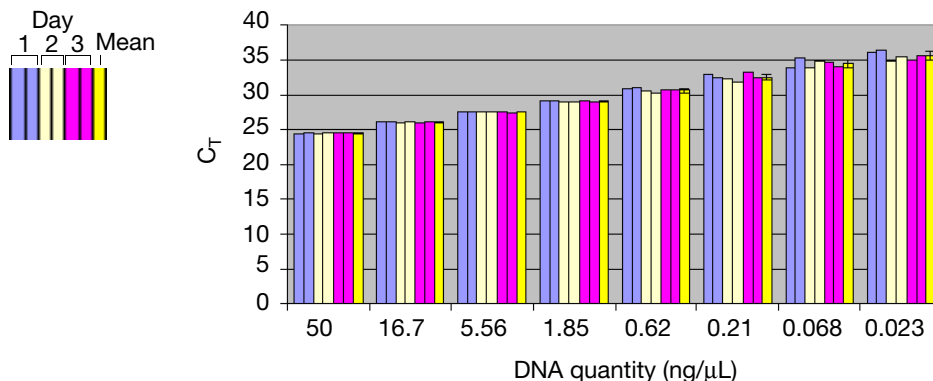


Figure 6-17 Precision: Quantifiler Y kit C_T results (7900HT SDS)

6.2.2 Mixture Studies (7900HT SDS)

An experiment was performed to demonstrate the specificity of the Quantifiler Human kit and the Quantifiler Y kit in analyzing mixtures of human genomic DNA from male and female sources. The mixture studies were designed to simulate circumstances in which a small component of male DNA must be discerned from a high background of female DNA.

Experiment

Purified genomic DNA from the Raji (male) and K-562 (female) cell lines were combined in ratios of 1:1, 1:4, 1:16, 1:64, 1:256 and 1:1024 (Raji:K-562). The male DNA was added at a constant level of 0.05 ng/μL in all samples, and the female DNA was present at amounts ranging from 0.05 ng/μL in the 1:1 sample to 50 ng/μL in the 1:1024 sample. The DNA amounts were based on the DNA concentrations provided by the suppliers and were not calibrated with the Quantifiler kits.

The mixtures were tested with the Quantifiler Human kit assay and the Quantifiler Y kit assay to determine the concentrations of total human genomic DNA (Quantifiler Human kit) and male DNA only (Quantifiler Y kit). For each sample, three replicate reactions were performed for each assay. Each assay used the same set of 8 human genomic DNA quantification standards run in duplicate reactions for each assay and both assays were run on the same reaction plate. The reaction plates were run on a 7900HT instrument.


Results The quantification results from using the Quantifiler Human kit varied from an average of 0.12 ng/μL for the 1:1 sample to 60 ng/μL for the 1:1024 sample, consistent with the increasing amounts of female DNA present.

The quantification results from using the Quantifiler Y kit varied from between 0.023 ng/μL to 0.058 ng/μL for all samples, regardless of the amount of female DNA present.

For the 1:1024 sample, the results showed a ratio of male DNA to total DNA of 1:1700. Differences between target concentrations and actual measurements were expected because the amounts of DNA added to the mixtures were based on the DNA concentrations provided by the suppliers and were not calibrated with the Quantifiler kits.

The results showed that the male DNA was detected and quantified accurately in all samples, regardless of the amount of female DNA present.

■ Quantifiler Y kit
■ Quantifiler Human kit

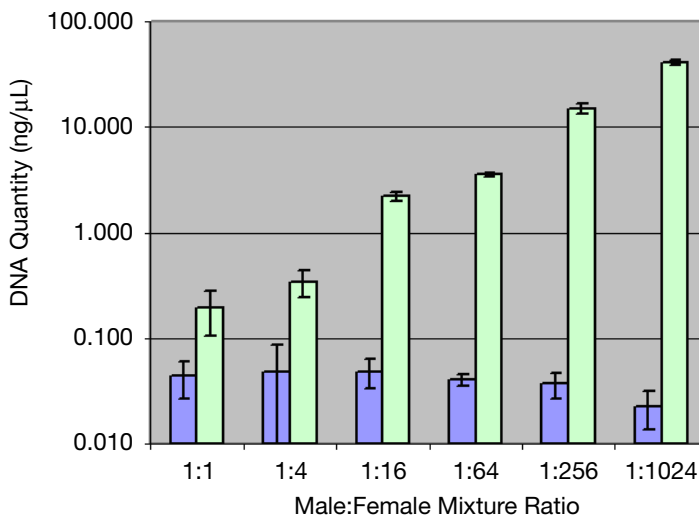

6.2.3 Comparisons with Other Methods (7900HT SDS)

Experiment Six human genomic DNA samples were obtained from different commercial sources.

Table 6-16 DNA samples tested with A_{260} and dye intercalation (7900HT SDS)

DNA	Sex	Extraction Source
007	Male	Blood
9948	Male	Cell line
Human genomic	Male	Blood
Raji-1	Male	Cell line
Raji-2	Male	Cell line
K-562	Female	Cell line

Using the concentrations provided by the supplier, the DNA samples were diluted to 2.0 ng/ μ L (A), 0.5 ng/ μ L (B), and 0.1 ng/ μ L (C).

Note: All dilutions were made in $T_{10}E_{0.1}$ Buffer with 20 μ g/mL glycogen added as a carrier and stabilizer.

All sample dilutions were quantified using the following methods:

- **A_{260}** – Because the concentrations of the dilutions extended below the detection limit of the spectrophotometer, absorbance at 260 nm was measured only for the highest dilution (2.0 ng/ μ L).

DNA concentration was calculated from the formula:

$$\text{Concentration } (\mu\text{g/mL}) = 50 \times A_{260}$$

The results calculated for the 2.0 ng/ μ L dilutions were then extrapolated for the higher dilutions (0.5 ng/ μ L and 0.1 ng/ μ L) using the known dilution factors.

- **Dye intercalation** – The microplate assay mode was used and the plate was read on a 7700 SDS. All of the sample dilutions were within the detection range of the assay. The assay was run using the λ bacteriophage DNA quantification standard supplied with the kit and a quantification standard based on Raji human genomic DNA. There were significant differences between the

standard curves from the λ DNA and Raji DNA. The results obtained from using the Raji DNA standard were used in these experiments because the Raji DNA was considered to be more similar to the DNA measured and because the results from using the Raji DNA standard were closer to the results obtained by the other methods.

- **Quantifiler kits** – DNA was quantified using the standard procedure. The Quantifiler Human DNA standard provided with the kits was used as recommended, with duplicate reactions for each of eight serial dilutions.

For each sample, the percent differences between Quantifiler kit results and results from the other two methods were calculated. The differences were expressed as a percentage of the reference method. For each method, the average percent differences from Quantifiler kit results were calculated. For comparisons with the Quantifiler Y kit, only results from male samples were used.

Results Table 6-17 shows the DNA concentrations calculated for all samples using the A_{260} method, the dye intercalation method, the Quantifiler Human kit, and the Quantifiler Y kit. It also shows the percent differences calculated for the comparisons.

Table 6-17 Comparison with A_{260} and dye intercalation (7900HT SDS)

Sample	A_{260} Result (ng/ μ L)	DI ^a Result (ng/ μ L)	Quantifiler Human Kit			Quantifiler Y Kit		
			Result (ng/ μ L)	% Diff. from A_{260}	% Diff. from DI	Result (ng/ μ L)	% Diff. from A_{260}	% Diff. from DI
007 A	2.74	2.502	2.094	23.6	16.3	3.547	29.4	41.8
007 B	0.685	0.756	1.007	47.0	33.2	0.950	38.6	25.7
007 C	0.137	0.176	0.272	98.8	54.6	0.220	60.3	24.6
9948 A	1.9	2.286	2.215	16.6	3.1	2.562	34.8	12.1
9948 B	0.475	0.496	0.677	42.5	36.4	0.634	33.5	27.8
9948 C	0.095	0.103	0.144	51.1	39.3	0.115	21.1	11.6
Human genomic A	2.2	2.270	2.887	31.2	27.2	1.792	18.6	21.1

Table 6-17 Comparison with A₂₆₀ and dye intercalation (7900HT SDS)

Sample	A ₂₆₀ Result (ng/ μ L)	DI ^a Result (ng/ μ L)	Quantifiler Human Kit			Quantifiler Y Kit		
			Result (ng/ μ L)	% Diff. from A ₂₆₀	% Diff. from DI	Result (ng/ μ L)	% Diff. from A ₂₆₀	% Diff. from DI
Human genomic B	0.55	0.584	0.805	46.3	37.9	0.379	31.0	35.0
Human genomic C	0.11	0.134	0.184	67.4	36.9	0.105	4.3	21.7
K-562 A	2.76	1.317	1.631	40.9	23.9	0.000	n.d.	n.d.
K-562 B	0.69	0.365	0.474	31.4	29.9	0.000	n.d.	n.d.
K-562 C	0.138	0.104	0.060	56.2	42.1	0.000	n.d.	n.d.
Raji-1 A	2	1.271	1.702	14.9	33.9	2.101	5.0	65.2
Raji-1 B	0.5	0.351	0.483	3.4	37.7	0.547	9.5	56.1
Raji-1 C	0.1	0.085	0.094	6.4	10.6	0.109	8.5	28.3
Raji-2 A	1.98	1.262	1.555	21.5	23.2	2.134	7.8	69.1
Raji-2 B	0.495	0.357	0.446	9.9	24.9	0.606	22.4	69.7
Raji-2 C	0.099	0.110	0.081	17.7	26.1	0.126	27.0	14.0

a. Dye intercalation method

The different methods produced similar quantification results.

Table 6-18 Average differences from A_{260} and dye intercalation (7900HT)

Method	Average Difference (%)	
	Quantifiler Human Kit	Quantifiler Y Kit
A_{260}	34.8	23.5
Dye intercalation	29.8	34.9

Section 6.3 Casework Sample Analysis

Case Type Studies

There is a recommended optimal DNA concentration range for using AmpF ℓ STR[®] PCR Amplification kits. The recommended amount of DNA input for the AmpF ℓ STR[®] Identifiler[®] PCR Amplification Kit is 0.5 to 1.25 ng human DNA (total per reaction), and for four-dye assays such as the AmpF ℓ STR[®] Profiler Plus[®] PCR Amplification Kit, 1.0 to 2.5 ng.

DNA quantification is specified as a requirement by the Scientific Working Group on DNA Analysis Methods (SWGDM) as a preliminary step to STR genotyping (Scientific Working Group on DNA Analysis Methods, 2000).

Experiment

A set of samples consisting of both non-casework and casework samples was tested. Of the sample set, 6 samples were non-casework, consisting primarily of blood sample extracts from single sources, and 22 were casework DNA extracts from fabric, clothing, or surface swabs. All DNA samples were prepared by organic extraction.

The DNA samples were quantified using the QuantiBlot[®] Human DNA Quantitation Kit (Applied Biosystems) and the Quantifiler Human kit performed on both the 7000 SDS and 7900HT SDS. The QuantiBlot kit was used in the chemiluminescent autoradiography mode. Tests with the Quantifiler kits for the 7000 SDS and 7900HT SDS were performed according to the standard procedure.

Using the results from the Quantifiler Human kit and the 7000 SDS, between 0.8 and 1.4 ng human genomic DNA was added to each Identifiler kit reaction, with many of the samples added at approximately 1.0 ng per reaction. Identifiler kit reactions were processed on the ABI PRISM[®] 3100 Genetic Analyzer and analyzed using GeneScan[®] Software v3.7.1 and Genotyper[®] Software v3.7, for use with the Windows NT[®] operating system. The STR profiles obtained from using the Identifiler kit were analyzed. Successful STR profiles produced complete profiles with peak heights between 200 and 4000 relative fluorescence units (RFU).

Results According to the results from the Quantifiler Human kit reactions run on the 7000 SDS, the range of DNA concentrations was 0.06 ng/μL to 2.61 ng/μL (Table 6-19).

Successful STR profiles were obtained for the 28 samples that were analyzed (Figure 6-18). These samples contained the minimum amount of DNA recommended for optimal Identifiler kit results (50 pg/μL in a 10-μL reaction). For some samples in the original set, the volume of DNA sample remaining after DNA quantification was insufficient to perform STR assays; these samples were not included in the data presented.


Table 6-19 Input for STR analysis of casework samples

STR ^a	Sample	Quantity (ng/μL)			Input Amount for Identifiler Kit (ng)
		QuantiBlot Kit	Quantifiler Human Kit and 7000 SDS	Quantifiler Human Kit and 7900HT SDS	
①	Non-casework	0.4	0.42	0.4	1.0
②	Non-casework	0.4	0.50	0.50	1.2
③	Non-casework	0.4	0.23	0.38	1.3
④	Non-casework	0.4	0.54	0.56	1.3
⑤	Non-casework	0.16	0.17	0.23	1.0
⑥	Non-casework	0.4	0.67	0.65	1.4
⑦	Positive control	n.d.	n.d.	n.d.	1.0
⑧	Negative control	n.d.	n.d.	n.d.	0.0
⑨	Cutting from shirt	0.4	0.78	0.88	1.1
⑩	Cutting from shirt	0.4	0.66	0.99	1.1
⑪	Cutting from fabric	0.06	0.093	0.11	1.2
⑫	Cutting from fabric	0.06	0.060	0.087	0.8

Table 6-19 Input for STR analysis of casework samples (*continued*)

STR ^a	Sample	Quantity (ng/μL)			Input Amount for Identifier Kit (ng)
		QuantiBlot Kit	Quantifiler Human Kit and 7000 SDS	Quantifiler Human Kit and 7900HT SDS	
13	Cutting from denim	0.16	0.10	0.13	1.3
14	Cutting from sock	0.04	0.11	0.15	1.1
15	Cutting from sweatshirt	1.2	2.61	3.75	1.4
16	Cutting from cotton	0.4	0.52	0.87	1.1
17	Cutting from sweatshirt	0.4	0.94	0.97	1.0
18	Cutting from cloth	0.4	0.31	0.56	1.1
19	Cutting from fabric	0.04	0.23	0.34	1.1
20	Cutting from leather	0.08	0.10	0.18	1.2
21	Cutting from carpet	0.4	0.76	0.95	1.3
22	Cutting from cloth	1.6	1.89	2.95	1.1
23	Cutting from shirt	1.2	2.29	3.10	1.2
24	Swab from hammer	0.6	0.47	0.58	1.1
25	Cutting from cloth	0.4	0.45	0.58	1.1
26	Cutting from fabric	0.08	0.16	0.18	1.3
27	Cutting from carpet	0.4	1.45	1.62	1.3
28	Cutting from cap	0.4	0.45	0.46	1.0
29	Cutting from shirt	1.2	2.29	3.10	1.3

a. See [Figure 6-18](#) for STR profiles.

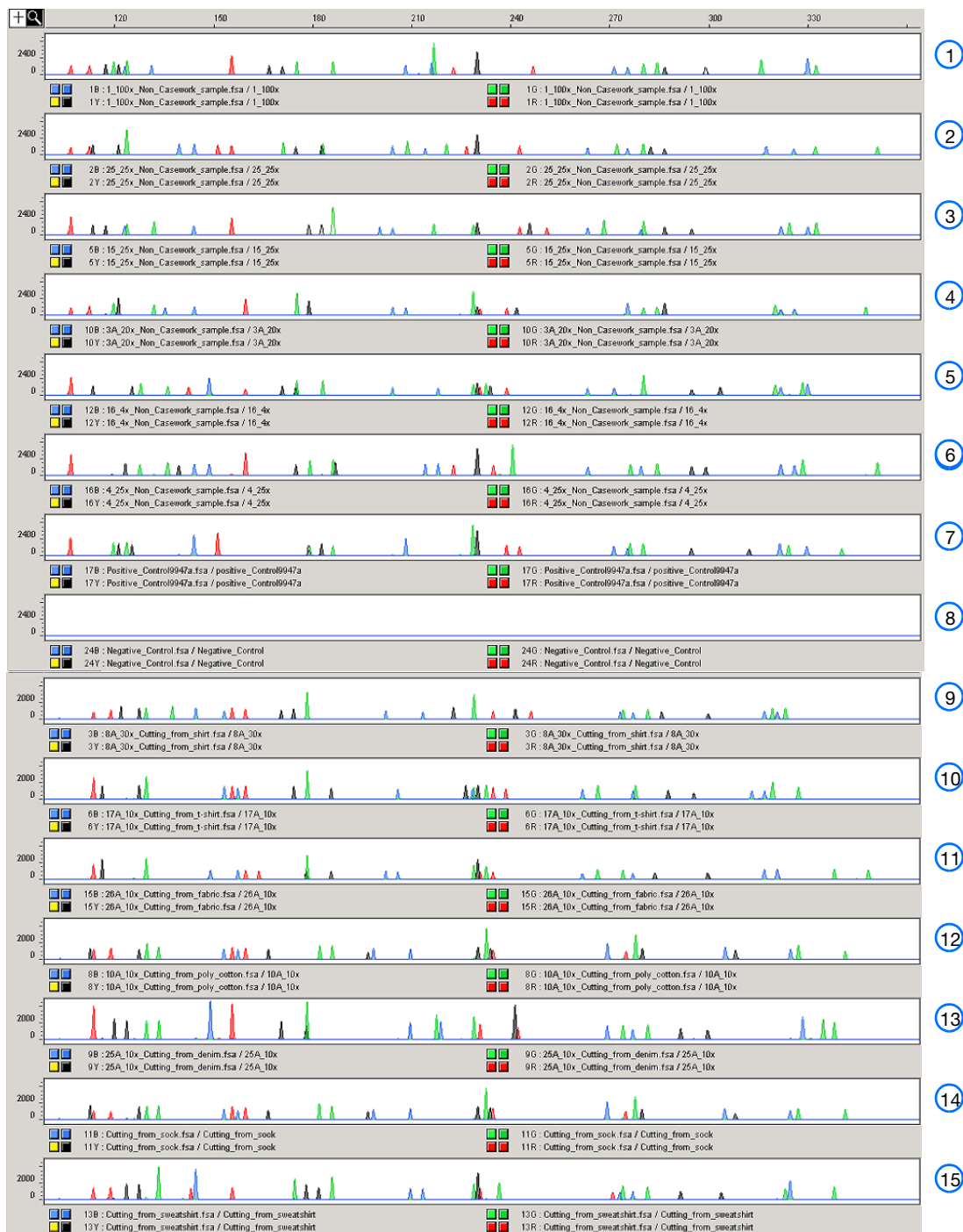

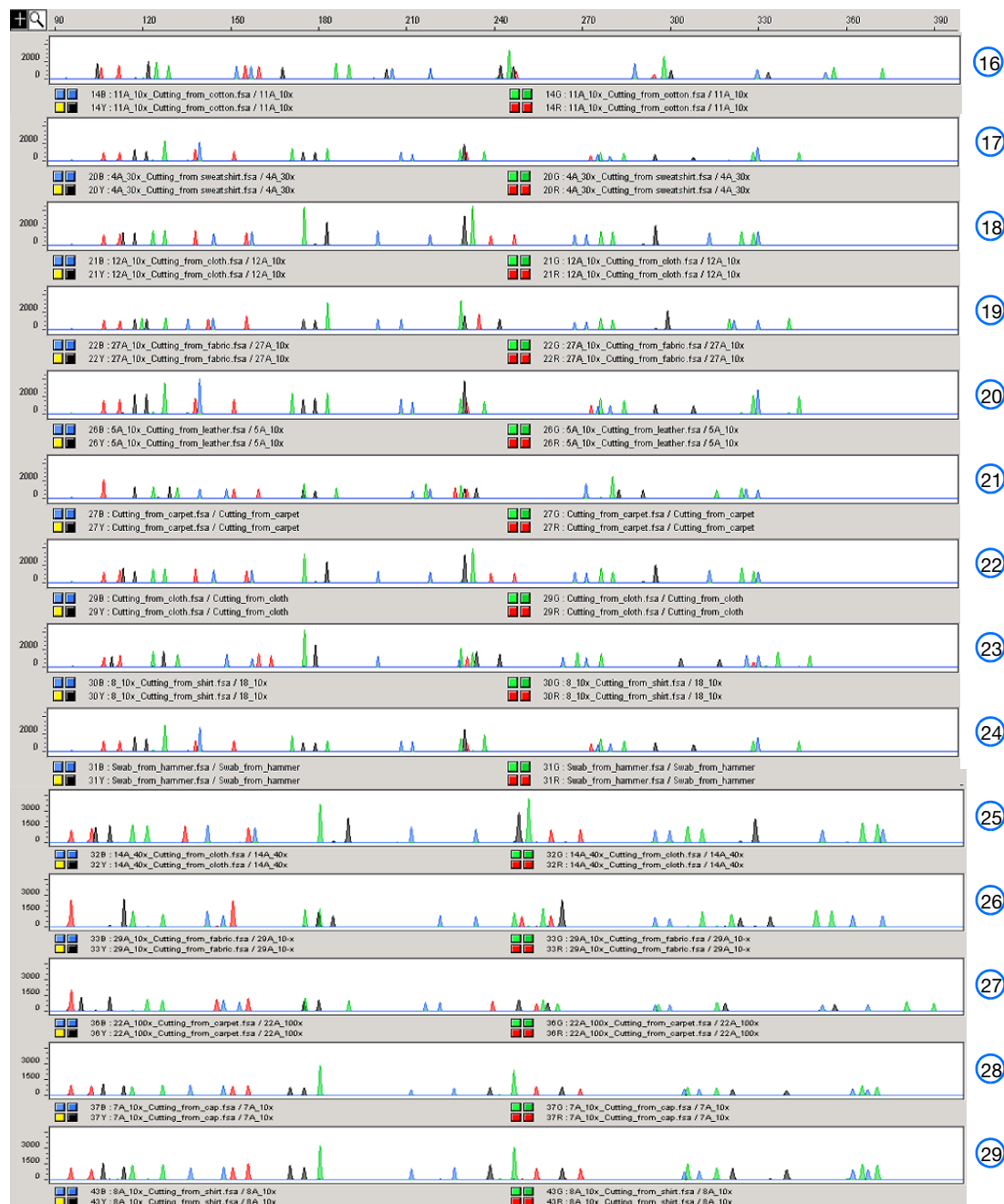


Figure 6-18 STR profiles of casework samples

Section 6.4 Applied Biosystems

7500 Real-Time PCR System Validation (SDS Software v1.2.3)

Overview The Quantifiler® Human DNA Quantification Kit and the Quantifiler® Y Human Male DNA Quantification Kit were tested (see the experiments listed below) using the Applied Biosystems 7500 Real-Time PCR System with SDS Software v1.2.3, running on the Windows® XP operating system. The results were then compared to the previously validated ABI PRISM® 7000 Sequence Detection System with SDS Software v1.0.

The experimental data generated demonstrate that the 7500 System (SDS Software v1.2.3):

- Provides accurate results when used with the Quantifiler kits for the analysis of genomic DNA samples.
- Produced results that are statistically similar to the results produced on the previously validated 7000 System (SDS Software v1.0).

Validation Experiments Performed

- Precision and Accuracy
- Reproducibility and Sensitivity
- Background
- Auto Baseline versus Manual analysis

6.4.1 Materials and Methods

6.4.1.1 Reagents

To minimize variables from hand pipetting and lot-to-lot reagent differences, the following set up procedures were used throughout the study:

- Eight serial dilutions were made with one lot of standard DNA provided with the Quantifiler kits (first dilution prepared with 500 µL DNA and 1,000 µL 10 mM Tris-HCl (pH 8.0) and 0.1 mM Na₂EDTA (T₁₀E_{0.1} buffer)).

- One manufactured lot of each kit was used for all validation studies:

Kit	Part Number	Lot Number
Quantifiler Human Kit	4343895	0501020
Quantifiler Y Human Male Kit	4343906	0501018

6.4.1.2 Instruments

Three 7500 systems (SDS Software v1.2.3) and three 7000 systems (SDS Software v1.0) were used for this study (six instruments total). Before the study, each instrument was calibrated by an Applied Biosystems service engineer (ROI calibration, background calibration, optical calibration, pure dye calibration, RNase P run).

The Biomek[®] FX Laboratory Automation Workstation was used to set up the real-time PCR reaction plates to minimize hand-pipetting variations:

- The PCR master mixes (PCR reagents with standard or sample DNA mixed together) were aliquoted into a 96-well plate (PCR master mix plate).
- Six empty 96-well plates and the PCR master mix plate were placed on the Biomek FX work surface.
- The Biomek FX aspirated 25 μ L from the PCR master mix plate, then slowly dispensed it into the corresponding well in an empty 96-well plate. The plates were sealed, spun down, then quickly loaded onto a 7500 or 7000 system. This process ensured timely and precise replication of real-time PCR plates for six instruments at a time.

6.4.2 Experimental Setup

Precision and Accuracy Testing

On each 96-well reaction plate, six sets of standard dilutions for each Quantifiler kit were set up for real-time PCR. [Figure 6-19](#) shows the experimental plate layout.

For each instrument, six replicate plates were run consecutively. The cycle threshold (C_T), R^2 , and slope values were compared statistically to determine precision and accuracy, which established 95% confidence intervals for each instrument type.

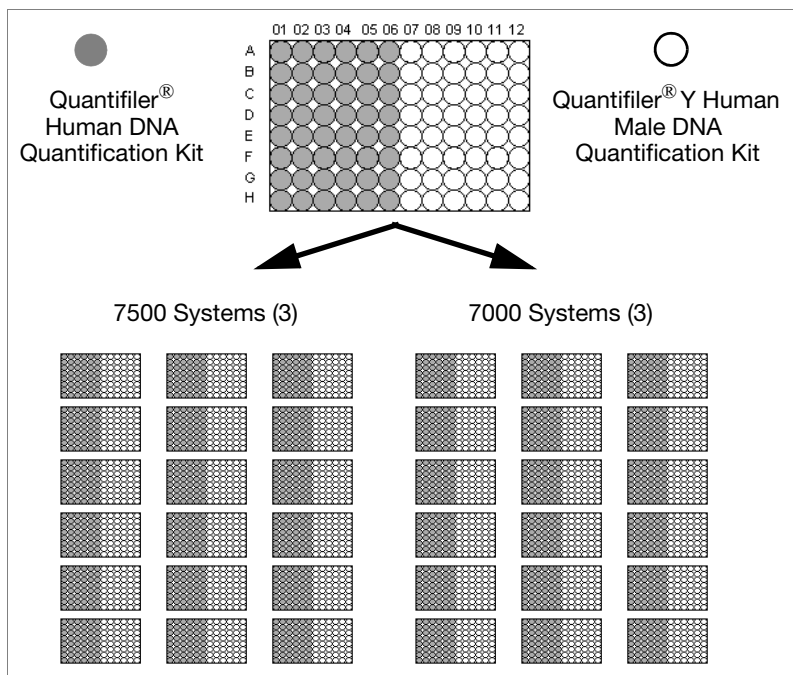


Figure 6-19 Plate layout – Precision and accuracy testing on the 7500 System (SDS Software v1.2.3) and 7000 System (SDS Software v1.0)

Reproducibility and Sensitivity Testing

On each 96-well reaction plate, the following were set up for real-time PCR:

- Standard dilution series (two replicates of each dilution point)
- Five replicate serial dilution sets of two sample DNAs (Raji and 9948B)

The experimental plate layout is shown in [Figure 6-20](#).

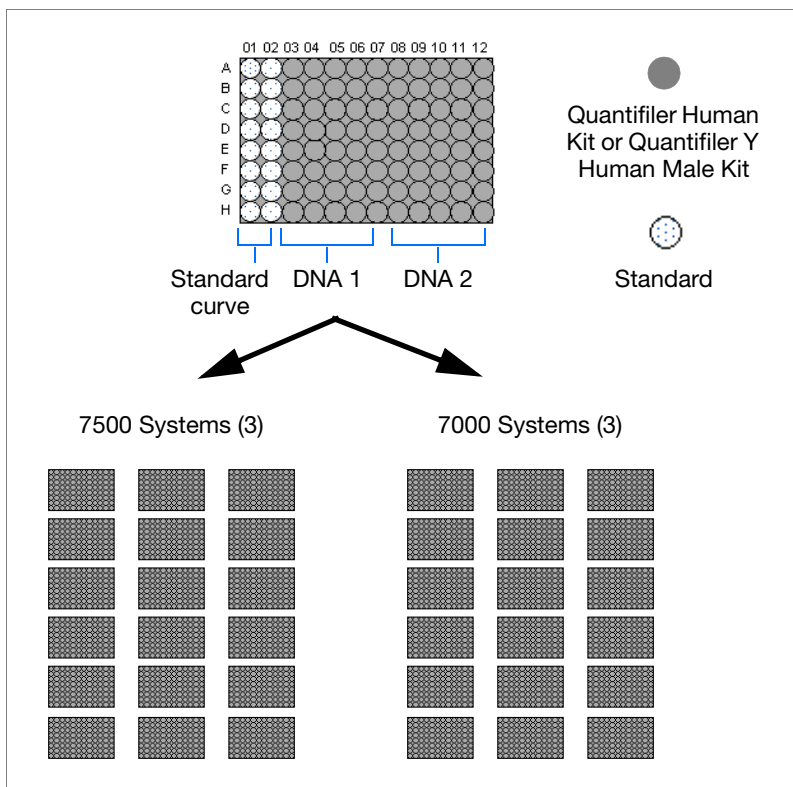


Figure 6-20 Plate Layout – Reproducibility and sensitivity testing on the 7500 System (identical plate layout for both kits)

On each instrument, six replicate plates were run consecutively with each Quantifiler kit (for a total of 18 plates on 7500 systems and 18 plates on 7000 systems).

To demonstrate reproducibility and sensitivity, the replicate DNA samples were quantitated, and the results were compared statistically between instrument types.

Background Testing Ninety-five no template controls (NTCs) and one positive control (the 50 ng/μL standard DNA dilution sample) were set up on a 96-well plate. One plate from each Quantifiler kit was run on each instrument (for a total of 12 plates).

6.4.3 Data Collection

The standard thermal cycling protocol (9600 Emulation mode) described in the [Chapter 3, “PCR Amplification,”](#) was used for all instrument runs.

6.4.4 Data Analysis

Initial Data Compiling and Analysis All runs were analyzed initially using Manual analysis mode, with the baseline set to 3 to 15 and the threshold set at 0.2.

Average values and standard deviations for C_T , slope, and R^2 were calculated for all replicate samples in a run.

For Auto-Baseline-to-Manual analysis comparisons, the run files from the 7500 System (SDS Software v1.2.3) were reanalyzed using Auto Baseline mode and a threshold of 0.2.

Statistical Data Analysis For statistical analysis, the Stat-Ease Design-Expert® Software was used for all ANOVA (analysis of variance) calculations. For paired t-Tests analysis, MicroSoft® Excel Analysis ToolPak software was used.

6.4.4.1 Precision and Accuracy

For the precision and accuracy tests between the two instrument types, the following values were determined:

- Average C_T
- Average Slope
- Average R^2
- 95% confidence intervals (CI) by ANOVA analysis

C_T Results Table 6-20 shows the average C_T values (95% CI) for the 7500 System (SDS Software v1.2.3) and the 7000 System (SDS Software v1.0) at each standard curve dilution.

Table 6-20 C_T Values (95% CI)

Standard Curve Dilution (ng/μL)	7500 System		7000 System	
	Average C _T Value (95% CI)	C _T Value Range (95% CI)	Average C _T Value (95% CI)	C _T Value Range (95% CI)
50	23.29	23.21 to 23.37	23.05	22.97 to 23.13
16.7	24.98	24.90 to 25.06	24.56	24.48 to 24.64
5.56	26.53	26.44 to 26.61	26.08	26.00 to 26.16
1.85	28.05	27.97 to 28.14	27.53	27.45 to 27.61
0.62	29.44	29.36 to 29.53	29.00	28.92 to 29.09
0.21	30.86	30.78 to 30.94	30.33	30.25 to 30.41
0.068	32.40	32.32 to 32.48	31.61	31.53 to 31.70
0.023	33.98	33.88 to 34.05	33.03	32.95 to 33.11

Statistically, the two instrument types resulted in significantly different C_T values ($p < 0.0001$) when compared with the ANOVA analysis. No significant difference in C_T values was observed when comparing results from instruments of the same type.

Slope Results Figure 6-21 shows the average slope values obtained for replicate standard curves run on each instrument. The slope values obtained for the 7500 System (SDS Software v1.2.3) are listed below and are within the ranges previously established on the 7000 System (SDS Software v1.0):

Kit	Slope	Established Slope Range
Quantifiler Human Kit	-2.93 to -3.18	-2.9 to -3.3
Quantifiler Y Human Male Kit	-3.05 to -3.36	-3.0 to -3.6

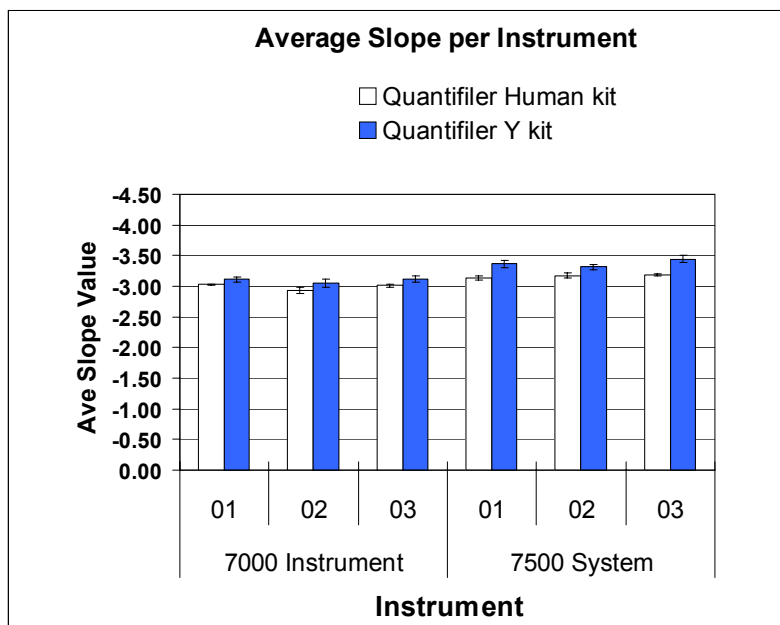


Figure 6-21 Average slope values – Replicate standard curves

R² Results Figure 6-22 shows the average R² values obtained for replicate standard curves on each instrument. All R² values were greater than 0.98 and are within the established range.

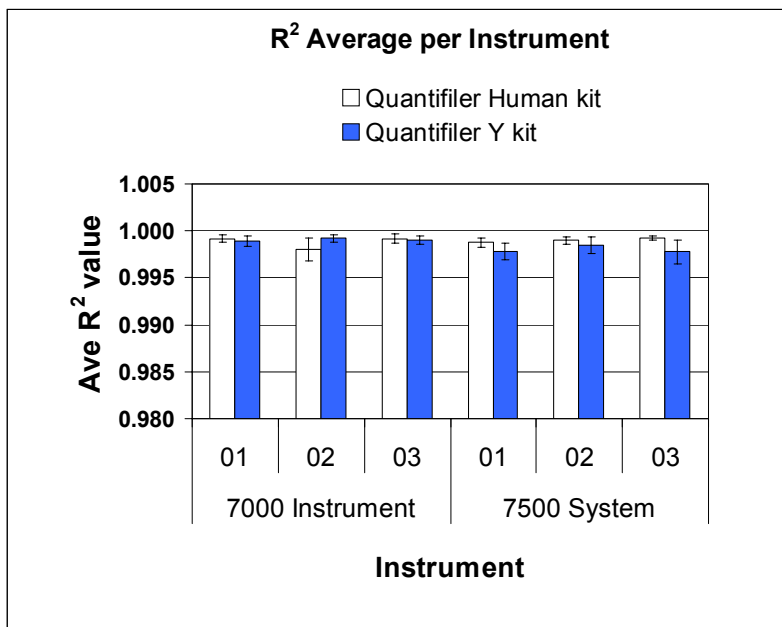


Figure 6-22 Average R² values – Replicate standard curves

6.4.4.2 Reproducibility and Sensitivity

Two sample DNAs were quantitated for this experiment. Eight 3-fold serial dilutions for each sample were run (five replicates per dilution, 40 wells per sample). The C_T values were generated in Manual analysis mode, then the quantities were calculated using the standard curve on each plate.

Figure 6-23 shows average C_T values (each point n = 90 replicates) across a set of four serial dilutions (2 ng/μL to 0.5 ng/μL) with the Quantifiler Human Kit and the corresponding quantitated concentrations for one DNA sample. Similar results were obtained for the second DNA sample and the Quantifiler Y Human Male Kit (data not shown).

As the data show, differences in C_T values do not affect calculated quantities (calculated quantities were normalized resulting in comparable concentrations on both instrument types).

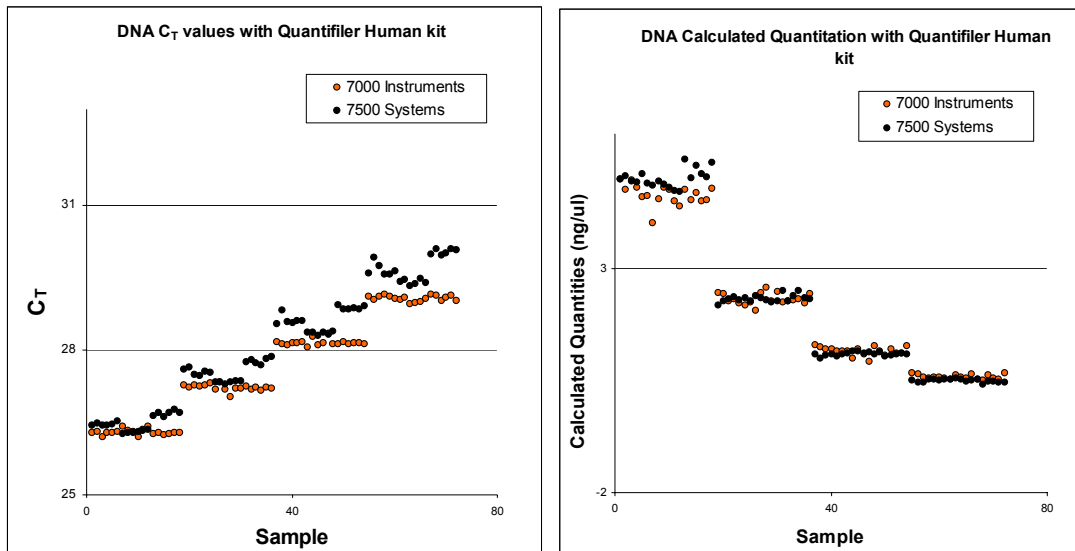


Figure 6-23 C_T values and quantitated concentrations – Quantifiler Human Kit (comparable data were obtained for the Quantifiler Y Human Male Kit)

Table 6-21 shows the average calculated quantities for each DNA sample obtained with the Quantifiler Human Kit. For sample concentrations between 2 ng/ μ L and 0.5 ng/ μ L, the percent difference between the quantitated values between instrument types did not exceed 16%. No statistically significant difference was observed for calculated quantities obtained using the Quantifiler Human Kit on the two instrument types.

Table 6-21 Average Calculated DNA Quantities – Quantifiler Human Kit

DNA Sample	7000 Avg Calculated Qty. (ng/ μ L)	7000 Std. Dev.	7500 Avg Calculated Qty. (ng/ μ L)	7500 Std. Dev.	Difference Between 7000 & 7500 Calculated Qty. (ng/ μ L)	% Difference of 7000 Qty. Value from 7500 Qty. Value
Raji	9.33	0.51	9.14	0.33	0.19	2.04
	4.58	0.15	4.24	0.12	0.34	7.72
	2.30	0.11	2.09	0.04	0.21	9.63
	1.16	0.05	1.07	0.03	0.09	8.01
	0.59	0.03	0.55	0.01	0.04	6.91
	0.27	0.01	0.26	0.01	0.01	3.43
	0.15	0.01	0.15	0.01	0.00	3.24
	0.08	0.00	0.07	0.00	0.01	8.04
9948	4.65	0.15	5.02	0.20	-0.37	7.58
	2.33	0.02	2.34	0.05	-0.01	0.36
	1.16	0.05	1.09	0.03	0.07	5.98
	0.59	0.02	0.50	0.03	0.08	15.52
	0.31	0.02	0.27	0.01	0.04	12.31
	0.17	0.01	0.15	0.01	0.02	10.80
	0.08	0.01	0.06	0.00	0.03	38.59
	0.05	0.01	0.04	0.00	0.01	18.14

Table 6-22 shows the average calculated quantities for each DNA sample obtained with the Quantifiler Y Human Male Kit. For sample concentrations of 2 ng/μL to 0.5 ng/μL, the percent difference between the quantitated values between instrument types did not exceed 18%. A minimal statistical difference was observed for calculated quantities obtained using the Quantifiler Y Human Male Kit on the two instrument types ($p = 0.0027$).


Table 6-22 Average Calculated DNA Quantities – Quantifiler Y Human Male Kit

DNA Sample	7000 Ave. Calculated Qty. (ng/μL)	7000 Std. Dev.	7500 Ave. Calculated Qty. (ng/μL)	7500 Std. Dev.	Difference Between 7000 & 7500 Calculated Qty. (ng/μL)	% Difference of 7000 Qty. Value from 7500 Qty. Value
Raji	9.12	0.40	9.09	0.07	0.03	0.34
	4.60	0.20	4.66	0.04	-0.06	1.29
	2.53	0.07	2.36	0.05	0.17	7.04
	1.29	0.09	1.19	0.03	0.10	8.12
	0.66	0.05	0.62	0.03	0.05	7.36
	0.33	0.02	0.30	0.02	0.02	7.89
	0.15	0.02	0.14	0.01	0.02	11.55
	0.070	0.02	0.057	0.01	0.01	19.85
9948	4.71	0.12	4.56	0.06	0.15	3.19
	2.43	0.14	2.30	0.06	0.12	5.14
	1.34	0.09	1.13	0.05	0.21	17.34
	0.68	0.03	0.62	0.03	0.06	9.93
	0.33	0.03	0.28	0.03	0.05	15.60
	0.18	0.01	0.14	0.01	0.04	24.87
	0.08	0.00	0.05	0.00	0.02	34.65
	0.04	0.00	0.03	0.00	0.01	38.29


6.4.4.3 Background

Figure 6-24 shows background amplification plots for 95 NTCs and one positive control for both kits (one plate each) run on the 7000 System (SDS Software v1.0). Figure 6-25 shows background amplification plots for the 7500 System (SDS Software v1.2.3).

On all instruments, the 95 NTC samples yielded negative results (all C_T values >40) with both Quantifiler kits.

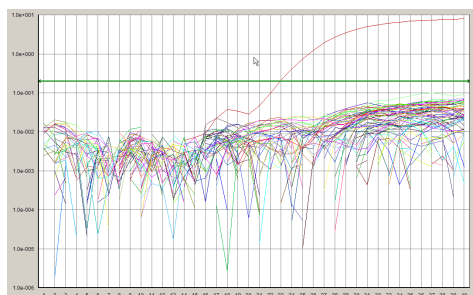


Quantifiler Human Kit

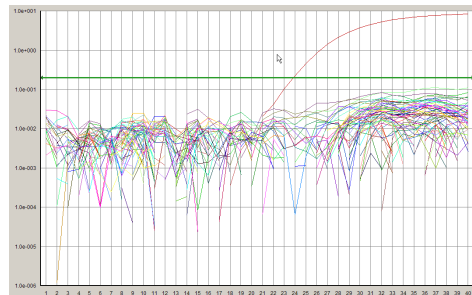


Quantifiler Y Human Male Kit

Figure 6-24 Background amplification plots – 7000 System (SDS Software v1.0)

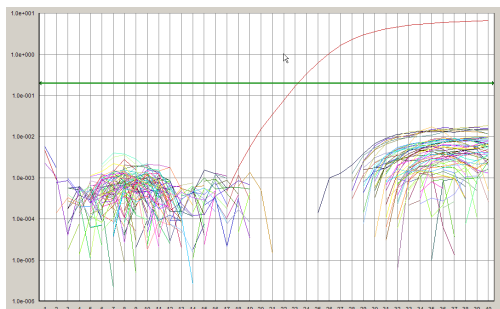

Quantifiler Human Kit

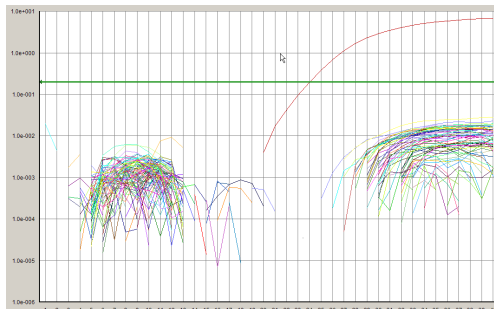

Quantifiler Y Human Male Kit

Figure 6-25 Background amplification plots – 7500 System (SDS Software v1.2.3)

6.4.4.4 Auto Baseline Analysis Versus Manual Analysis

C_T Precision and Accuracy

For Auto-Baseline-to-Manual analysis comparisons:

- The SDS software v1.2.3 data from the experiments described on the previous pages were reanalyzed in Auto Baseline mode (default threshold 0.2).
- The C_T values were compared to each other.

Figure 6-26 shows the C_T values obtained using the Auto Baseline and Manual analysis modes with the Quantifiler Human Kit. Similar data were obtained for the Quantifiler Y Human Male Kit.

No statistically significant differences were observed for C_T values generated using the Auto Baseline and Manual analysis modes with either Quantifiler kit.

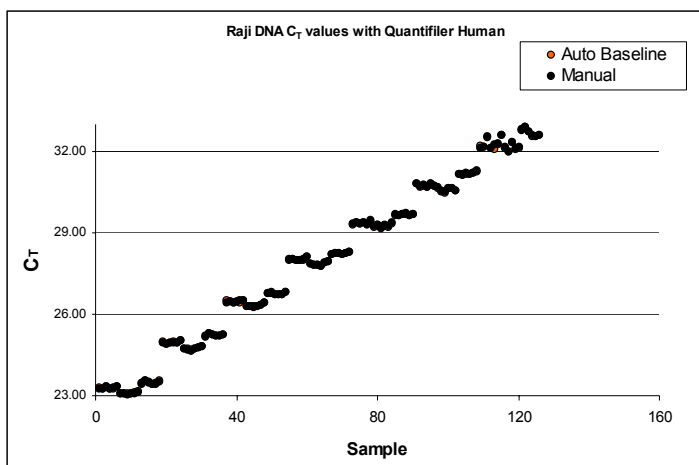


Figure 6-26 Comparison of C_T values between Auto Baseline and Manual analysis modes

C_T Reproducibility and Sensitivity

Figure 6-27 shows the C_T values and calculated quantities obtained using the Auto Baseline and Manual analysis modes with the Quantifiler Human Kit. Similar data were obtained for the Quantifiler Y Human Male Kit.

No statistically significant differences were observed for C_T values and calculated quantities derived using the Auto Baseline and Manual analysis modes with either Quantifiler kit.

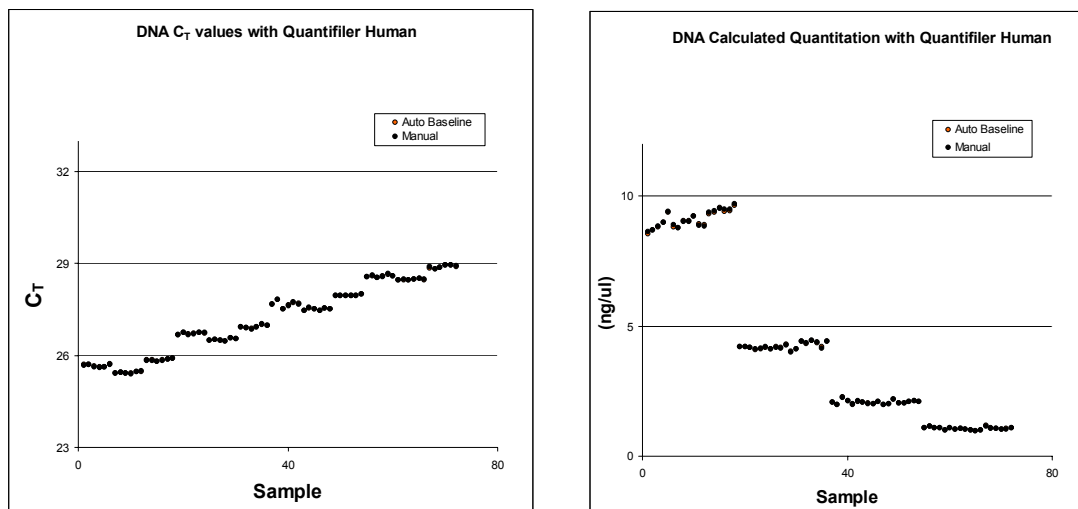


Figure 6-27 Comparison of C_T values and the corresponding calculated quantities – Auto Baseline and Manual analysis modes – Quantifiler Human Kit

6.4.5 Discussion

6.4.5.1 Precision and Accuracy

7500 System Comparison: No statistically significant differences were observed in C_T , slope, and R^2 values between replicate samples run on the 7500 System (SDS Software v1.2.3) using both Quantifiler kits.

7500-to-7000 System Comparison: Statistically significant differences in C_T , slope, and R^2 values were observed in samples run on the 7500 System (SDS Software v1.2.3) versus the 7000 System (SDS Software v1.0) using both Quantifiler kits.

However, the data obtained from both instrument types are within the previously established parameter ranges published in the *Quantifiler*® *User's Manual*, Chapter 5, Table 5-1.

6.4.5.2 Reproducibility and Sensitivity

Sensitivity: Similar C_T values and calculated DNA quantities were obtained at each of the standard curve concentrations, demonstrating similar sensitivity results between the 7000 System (SDS Software v1.0) and 7500 System (SDS Software v1.2.3).

Calculated Quantities: Data obtained using the Quantifiler Human Kit showed no statistically significant difference when the calculated quantities obtained from the 7000 and 7500 systems were compared ($p = 0.22$, with 95% confidence). However, minimally significant differences were observed between the two instrument types for calculated quantities using the Quantifiler Y Human Male Kit.

To further explore the extent of the difference between the two instrument types, the percent differences between the calculated quantities within the concentration range of 2 ng/ μ L to 0.5 ng/ μ L were determined. This range was selected because it represents the optimal input range for most STR kits. In this range, there was, at most, an 18% concentration difference between calculated quantities using the 7000 and the 7500 systems. The impact of the slight differences in calculated quantities should have minimal effect on results of STR analysis. However, laboratories should perform the appropriate studies to verify optimal input amounts for amplification.

6.4.5.3 Auto Baseline Analysis Versus Manual Analysis

No statistically significant difference was observed for C_T values and calculated quantities derived using the Auto Baseline and Manual analysis modes on the 7500 System (SDS Software v1.2.3).

6.4.6 Conclusion

This validation study demonstrates that the Applied Biosystems 7500 Real-Time PCR System with SDS Software v1.2.3 is a robust, reliable, and reproducible system for performing DNA quantification using the Quantifiler kits.

When statistically comparing 7500 System (SDS Software v1.2.3) results (C_T , slope, and R^2 values) to results obtained using previously validated ABI PRISM[®] 7000 Sequence Detection System with SDS Software v1.0:

- Differences in calculated quantities are minimal (Quantifiler Y Human Male Kit) or insignificant (Quantifiler Human Kit) for unknown samples using the 7500 and 7000 systems.
- The differences observed should have little effect on resulting STR amplification based on calculated DNA quantities.
- No significant difference is observed between C_T values and calculated quantities derived by using Auto Baseline and Manual analysis modes.

Section 6.5 ABI PRISM® 7000 Sequence Detection System Validation (SDS Software v1.2.3)

Overview The Quantifiler® Human DNA Quantification Kit and Quantifiler® Y Human Male DNA Quantification Kit were tested (see the experiments listed below) using the ABI PRISM® 7000 Sequence Detection System with SDS Software v1.2.3, running on the Windows® 2000 operating system., then compared to the previously validated ABI PRISM® 7000 Sequence Detection System with SDS Software v1.0.

The experimental data generated demonstrate that the 7000 System (SDS Software v1.2.3):

- Provides accurate results when used with the Quantifiler kits for the analysis of genomic DNA samples.
- Produced results that are similar to the results produced on the previously validated 7000 System (SDS Software v1.0).

Validation Experiments Performed

- Precision and Accuracy
- Reproducibility and Sensitivity
- Background
- Auto Baseline versus Manual analysis

6.5.1 Materials and Methods

6.5.1.1 Reagents

To minimize variables from hand-pipetting and lot-to-lot reagent differences, the following set-up procedures were used throughout the study:

- Eight serial dilutions were made with one lot of standard DNA provided with the Quantifiler kits (first dilution prepared with 500 µL DNA and 1,000 µL 10mM Tris-HCl (pH 8.0) and 0.1 mM Na₂EDTA (T₁₀E_{0.1} buffer)).

- One manufactured lot of each kit was used for all validation studies:

Kit	Part Number	Lot Number
Quantifiler Human Kit	4343895	0501022
Quantifiler Y Human Male Kit	4343906	0501020

6.5.1.2 Instruments

One ABI PRISM® 7000 Sequence Detection System was used for this study under the following conditions:

- All experiments were run initially using SDS Software v1.0.
- The 7000 system computer was upgraded to SDS Software v1.2.3.
- The 7000 System (SDS Software v1.2.3) was calibrated by an Applied Biosystems service engineer (background calibration, pure dye calibration, RNase P run).
- For the following experiments, v1.0 data was reanalyzed using SDS Software v1.2.3:
 - Precision and Accuracy
 - Reproducibility and Sensitivity
 - Background
- For Auto Baseline versus Manual analysis experiments, new data were collected using SDS Software v1.2.3, analyzed in Auto Baseline mode, then reanalyzed in Manual mode.

6.5.2 Experimental Setup

Precision and Accuracy Testing

On each 96-well reaction plate, six sets of standard dilutions for each Quantifiler kit were set up for real-time PCR. The experimental plate layout is shown in [Figure 6-28](#).

Three replicate plates were run consecutively. The C_T , slope, and R^2 values were compared to determine precision and accuracy.

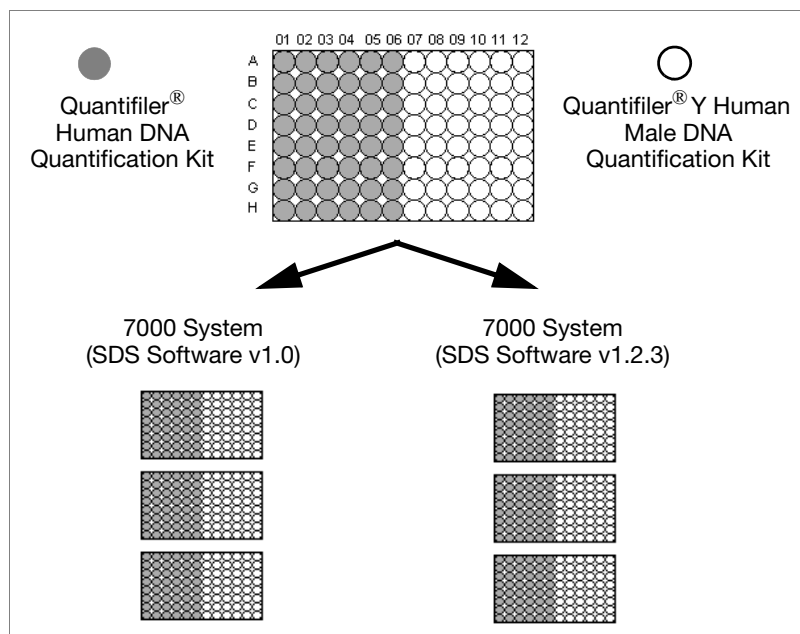


Figure 6-28 Plate Layout – Precision and accuracy testing on the 7000 System

Reproducibility Sensitivity, and Background Testing

On each 96-well reaction plate, the following were set up for real-time PCR:

- Standard dilution series (two replicates of each dilution point)
- Four replicate serial dilution sets of two sample DNAs (007 and 9948B)
- Sixteen no template controls (NTCs), which served as background samples

Figure 6-29 shows the experimental plate layout.

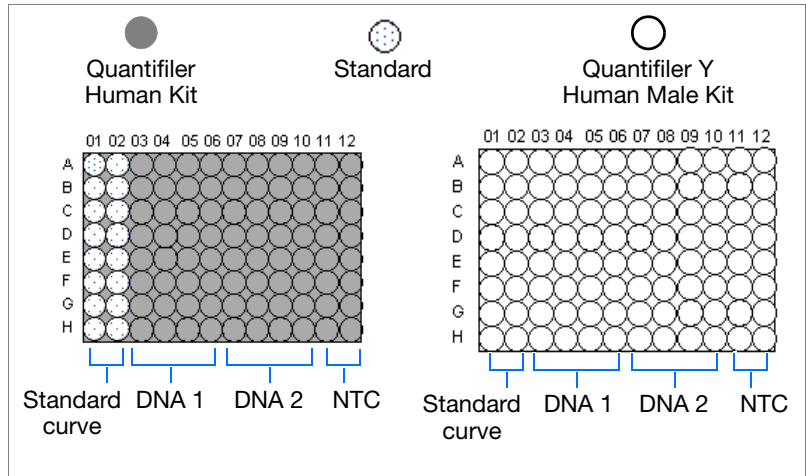


Figure 6-29 Plate Layout – Reproducibility and sensitivity experiments – 7000 Systems

One plate was run with each type of Quantifiler kit.

To demonstrate reproducibility and sensitivity, the replicate DNA samples were quantitated and the results were compared between each software version.

6.5.3 Data Collection

The standard thermal cycling protocol (9600 Emulation mode) described in the [Chapter 3, “PCR Amplification,”](#) was used for both studies.

6.5.4 Data Analysis

Initial Data Compiling and Analysis

All runs were analyzed initially using Manual analysis mode, with the baseline set to 3 to 15 and the threshold set at 0.2.

Average values and standard deviations for C_T , slope, and R^2 were calculated for all replicate samples in a run.

The instrument was then upgraded to SDS Software v1.2.3, then the same run files were reanalyzed and exported with the same analysis settings.

For Manual-to-Auto-Baseline analysis comparisons, the run files from the 7000 System (SDS Software v1.2.3) were reanalyzed using the Auto Baseline mode and a threshold of 0.2.

6.5.4.1 Precision and Accuracy

For the precision and accuracy tests between the two software versions, the average C_T , average slope, and average R^2 values were determined.

C_T Results Figures 6-30 to 6-32 show C_T values obtained using the SDS Software v1.0 and v1.2.3. The data consistently show that SDS Software v1.2.3 yields lower C_T values (2% difference).

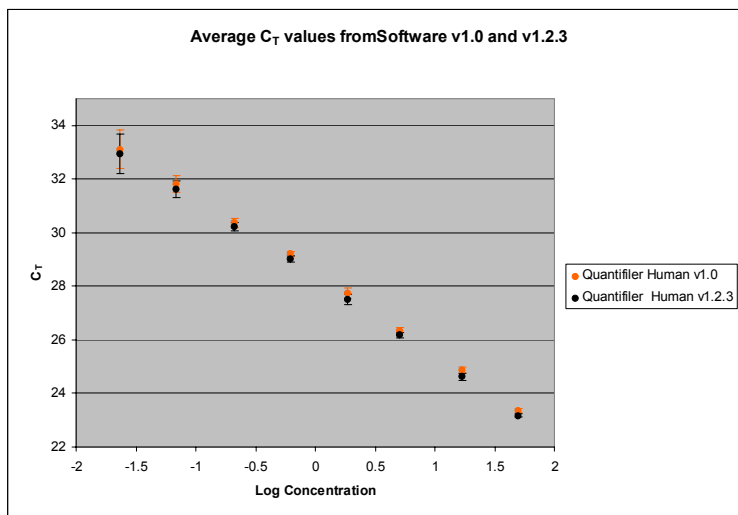


Figure 6-30 Average C_T values – Quantifiler Human Kit – SDS Software v1.0 and v1.2.3 (error bars indicate standard deviations)

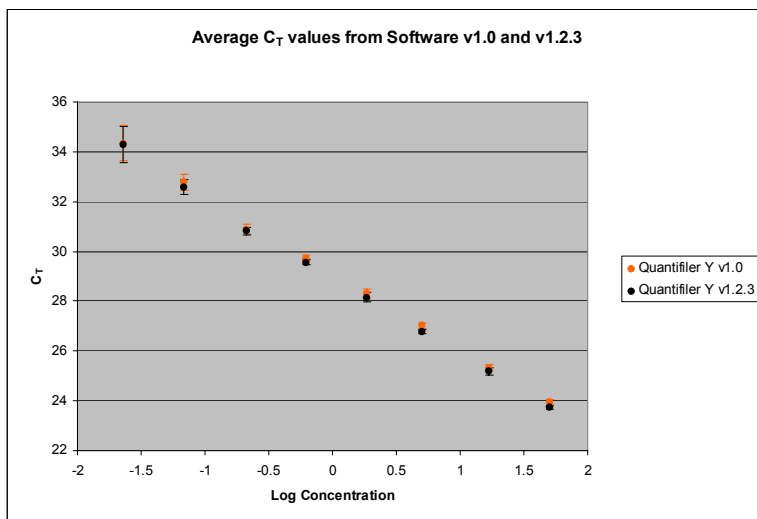


Figure 6-31 Average C_T values – Quantifiler Y Human Male Kit – SDS Software v1.0 and v1.2.3 (error bars indicate standard deviations)

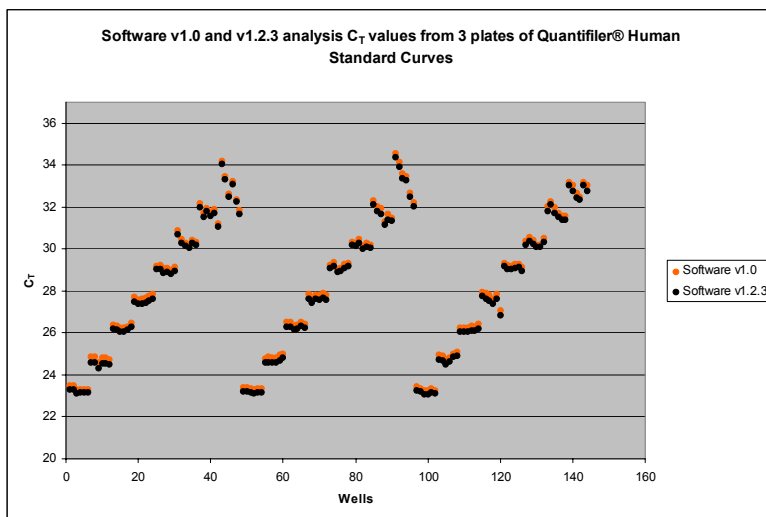


Figure 6-32 C_T Values per Sample – v1.0 compared to v1.2.3 – Quantifiler Human Kit

Slope Results Figure 6-33 shows the average slope values obtained using the SDS software v1.2.3 compared to v1.0. The slope values obtained for the 7000 System (SDS Software v1.2.3) are within the established ranges.

Kit	Slope	Established Slope Range
Quantifiler Human Kit	-2.90 to -2.97	-2.9 to -3.3
Quantifiler Y Human Male Kit	-3.0 to -3.09	-3.0 to -3.6

A 1% slope difference is observed between the v1.2.3 and v1.0 software.



Figure 6-33 Average slope values – SDS Software v1.0 and v1.2.3

R² Results Figure 6-34 shows that SDS software v1.2.3 yields data that are within the acceptable range of R² values: 0.98 to 1 for both kits (<0.5% difference).

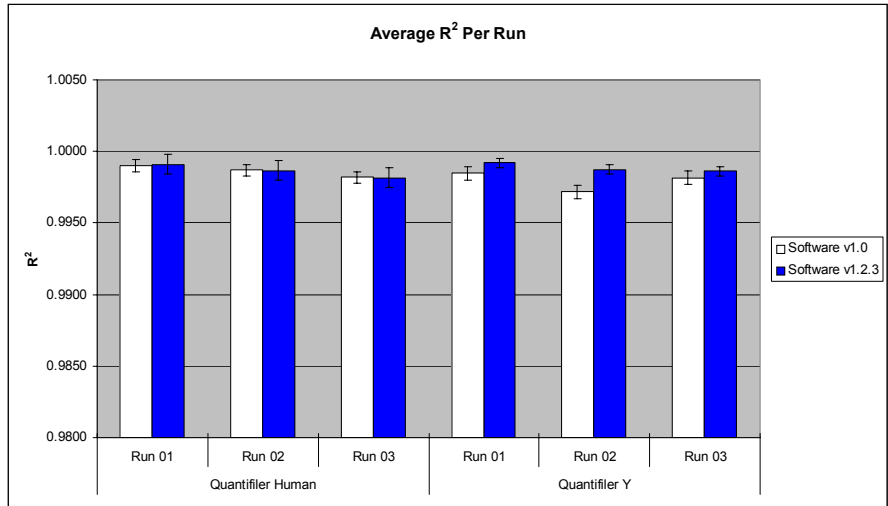


Figure 6-34 Average R² values – Quantifiler Human Kit and Quantifiler Y Human Male Kit – SDS Software v1.0 and v1.2.3.

6.5.4.2 Reproducibility and Sensitivity

Two sample DNAs were quantitated for this experiment. Eight 2-fold serial dilutions for each sample were run (four replicates per dilution, 32 wells per sample). The C_T values were generated in Manual analysis mode, then the quantities were calculated using the standard curve on each plate.

Figure 6-35 shows the C_T values for 007 and 9948B across a set of eight serial dilutions (~30 ng/μL to 0.1 ng/μL) with the Quantifiler Human Kit and the corresponding quantitated concentrations.

As the data show, differences in C_T values do not affect calculated quantities (calculated quantities were normalized, resulting in comparable concentrations from results generated with both software versions.)

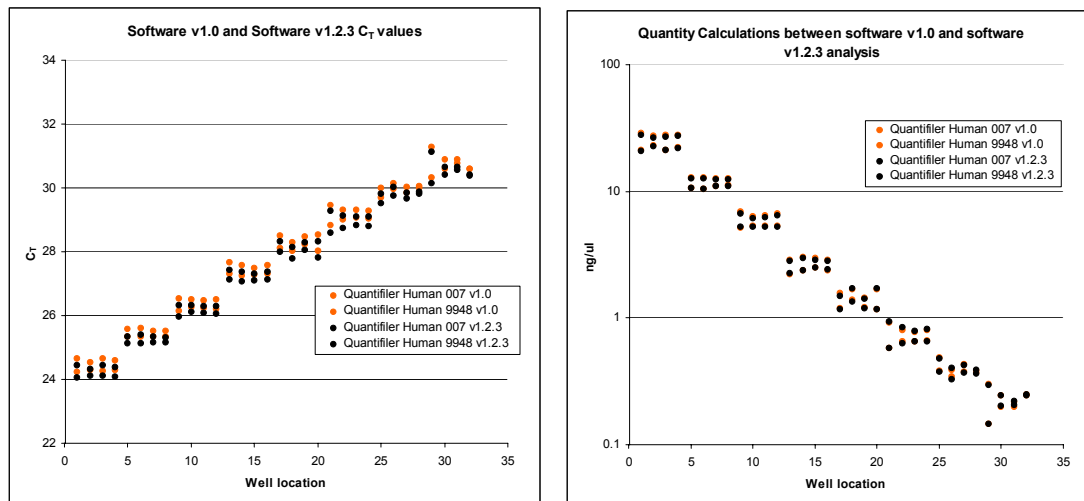


Figure 6-35 Average C_T values and quantitated DNA concentrations – 007 and 9948B – Quantifiler Human Kit

Figure 6-36 shows C_T results for the Quantifiler Y Human Male Kit that differ slightly between the v1.0 analysis and the v1.2.3 analysis. However, differences in C_T values do not affect calculated quantities (calculated quantities were normalized resulting in comparable concentrations from results generated with both software versions.)

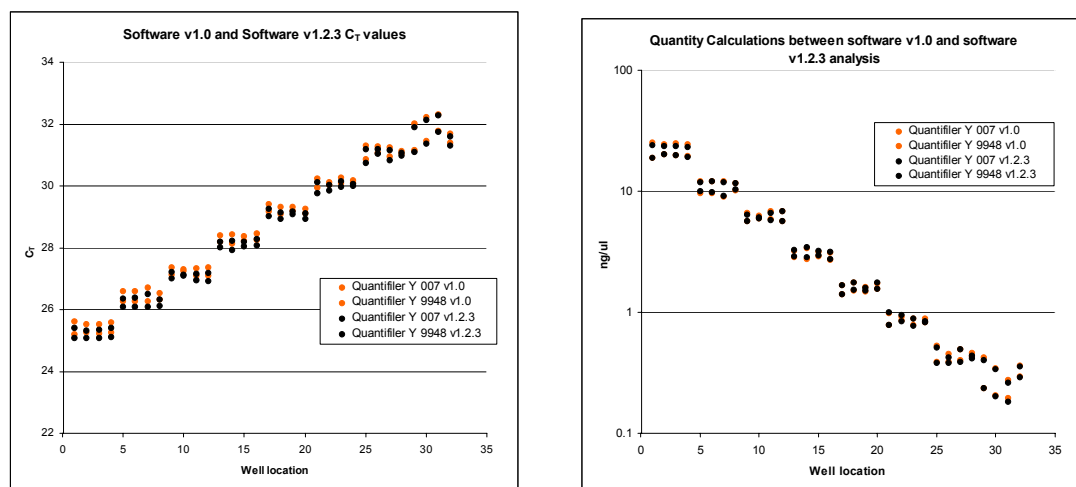


Figure 6-36 Average C_T values and quantitated DNA concentrations – 007 and 9948B – Quantifiler Y Human Male Kit

Figure 6-37 shows that there was a $\leq 6\%$ quantity difference between results obtained with v1.0 and v1.2.3 software.


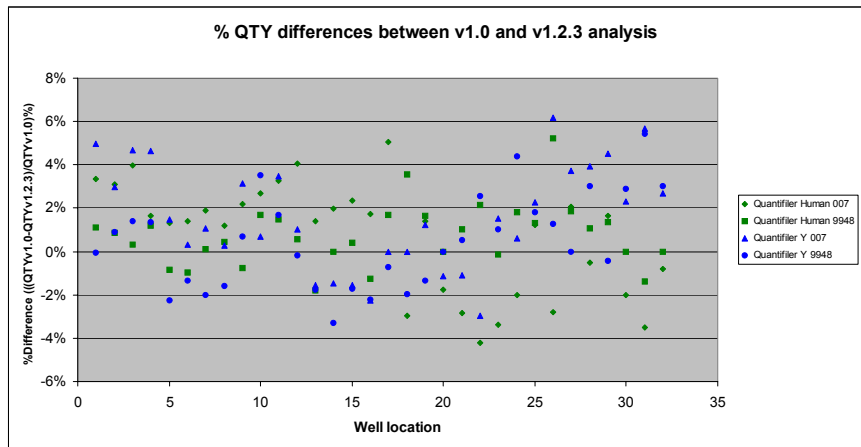



Figure 6-37 Percent DNA quantity differences – SDS Software v1.0 and v1.2.3

6.5.4.3 Background

Figure 6-38 shows the background results for 16 NTCs and one positive control for both kits run on the 7000 System (SDS Software v1.0). One out of 16 NTCs for the Quantifiler Human Kit resulted in a < 40 C_T result (36.81 C_T). Remaining NTCs resulted in > 40 C_T values (negative results).

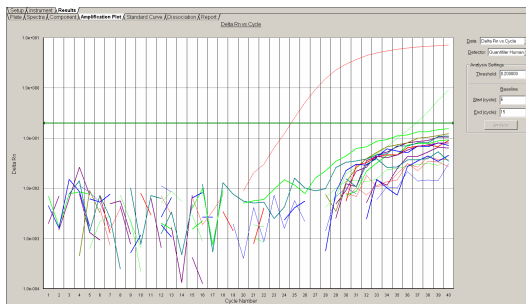


Quantifiler Human Kit

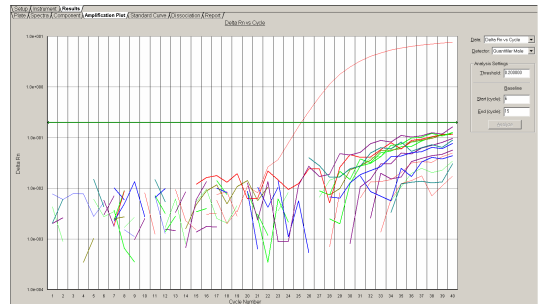


Quantifiler Y Human Male Kit

Figure 6-38 Background amplification plots – 7000 System (SDS Software v1.0)

Figure 6-39 shows the background results for 16 NTCs and one positive control for both kits reanalyzed on the 7000 System (SDS Software v1.2.3). One out of 16 NTCs for the Quantifiler Human Kit resulted in a <40 CT value ($38.26 C_T$). Overall, the NTC results do not change when analyzed with version 1.2.3.

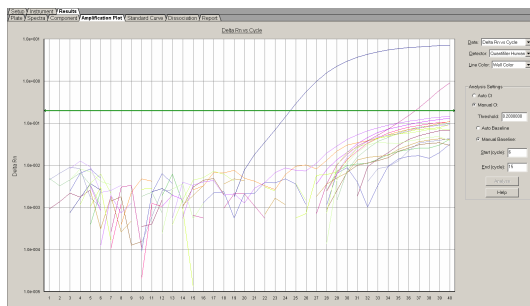

Quantifiler Human Kit

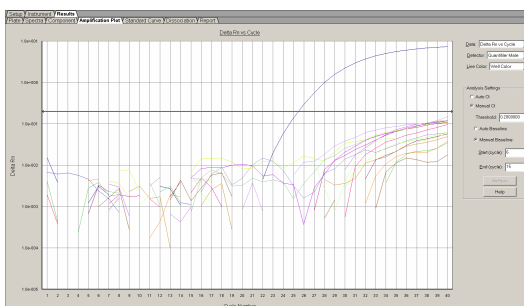

Quantifiler Y Human Male Kit

Figure 6-39 Background amplification plots – 7000 System (SDS Software v1.2.3)

6.5.4.4 Auto Baseline Analysis Versus Manual Analysis

C_T Precision and Accuracy

For Manual-to-Auto-Baseline analysis comparisons:

- Data from initial runs were collected with SDS Software v 1.2.3 and analyzed in Manual analysis mode, then reanalyzed in Auto Baseline analysis mode (default threshold 0.2).
- The C_T values were compared to each other.

Figures 6-40 and 6-41 show the average C_T values between Auto Baseline analysis and Manual analysis. There is a $<2\%$ difference between the two analysis methods for both kits.

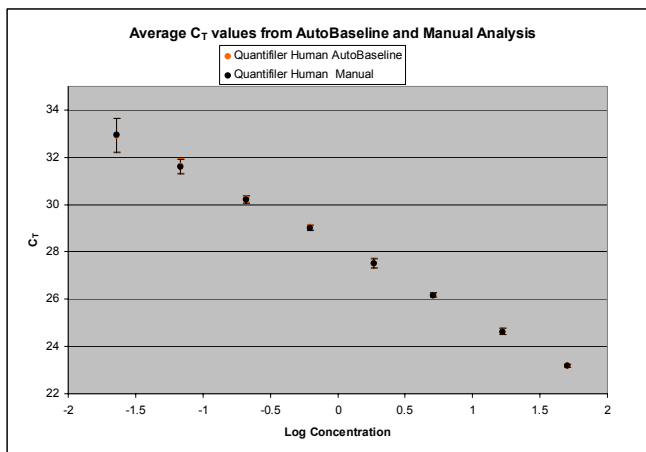


Figure 6-40 Comparison of C_T values between Auto Baseline and Manual analysis – Quantifiler Human Kit (error bars indicate standard deviations)

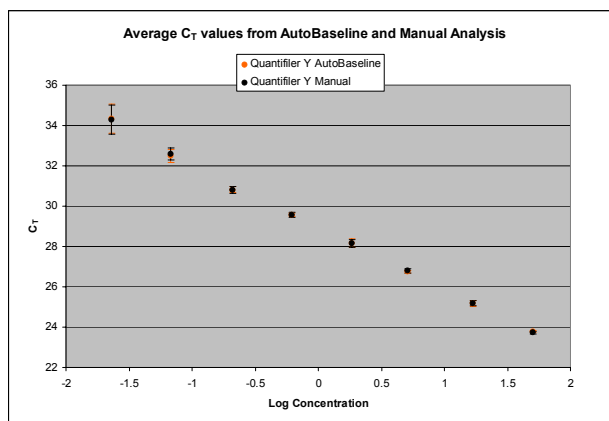


Figure 6-41 Comparison of C_T values between Auto Baseline and Manual analysis – Quantifiler Y Human Male Kit (error bars indicate standard deviations)

Reproducibility and Sensitivity

C_T

Figure 6-42 shows the C_T values obtained using the Auto Baseline and Manual analysis modes with the Quantifiler Human Kit.

No significant differences were observed for C_T values generated using the Auto Baseline and Manual analysis modes with either Quantifiler kit.

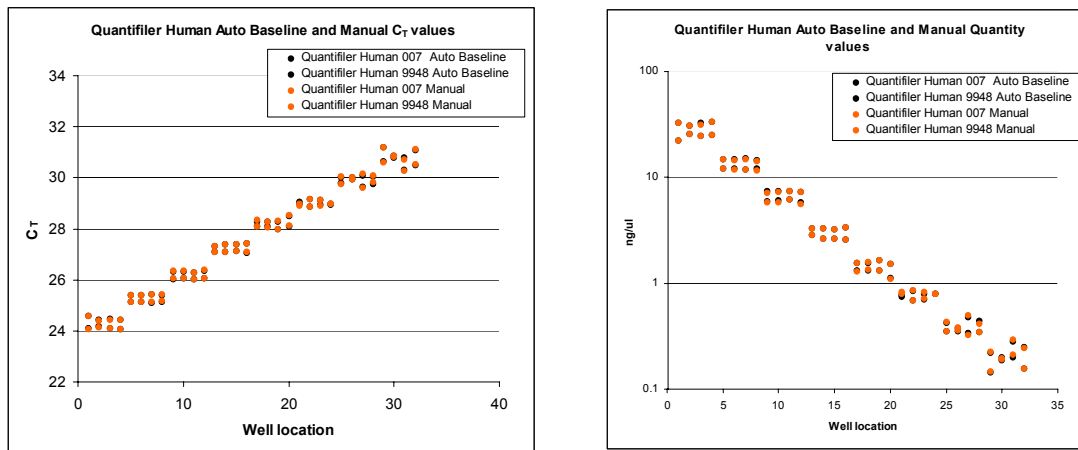


Figure 6-42 Average C_T values and average calculated quantities for 9948 and 007 – Quantifiler Human Kit (~30 ng/μL to 0.1 ng/μL)

Figure 6-43 shows the C_T values obtained using the Auto Baseline and Manual analysis modes with the Quantifiler Human Kit.

No significant differences were observed for C_T values generated using the Auto Baseline and Manual analysis modes with either Quantifiler kit. Auto Baseline C_T values overlap the manual C_T values. The corresponding quantities also overlap.

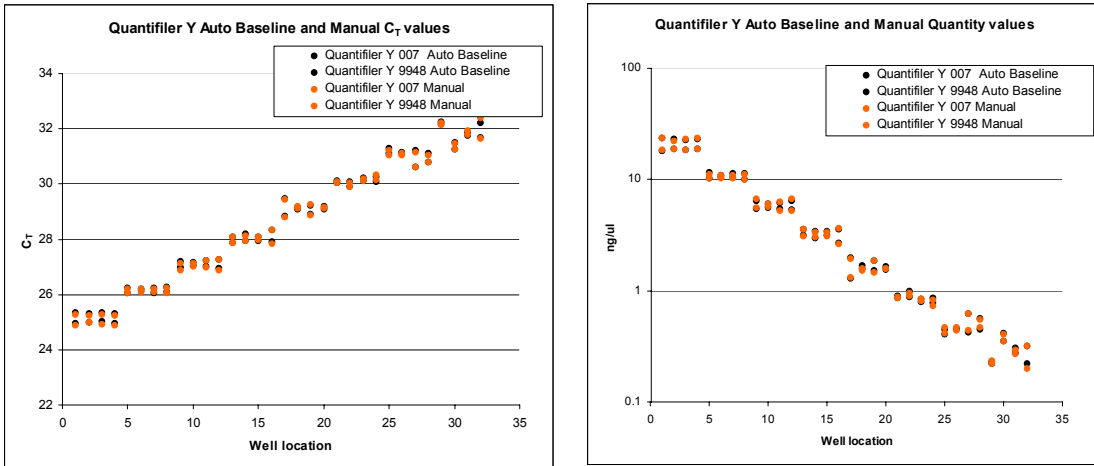


Figure 6-43 Average C_T values and average calculated quantities for 9948 and 007 – Quantifiler Y Human Male Kit

6.5.5 Discussion

6.5.5.1 Precision and Accuracy

The results from SDS Software v1.0 and v1.2.3 on a 7000 System slightly differ in C_T value (2% difference), slope (1%), and R^2 (<0.5%) for both Quantifiler kits. All v1.0 data and v1.2.3 data are within the *Quantifiler® User's Manual* published parameter ranges.

6.5.5.2 Reproducibility and Sensitivity

For both Quantifiler kits, there was a maximum difference of 6% when the calculated quantities using v1.0 and v1.2.3 were compared. Such minor differences in calculated quantities should not affect the ability to obtain interpretable STR profiles using the optimal input amount determined by individual laboratories during validation of the Quantifiler kits.

6.5.5.3 Manual Analysis Versus Auto Baseline Analysis

C_T values and their corresponding calculated quantities showed a maximum 8% difference between Auto Baseline and Manual analysis modes on the 7000 System (SDS Software v1.2.3). However, the differences observed should have little effect on resulting STR amplification based on calculated DNA quantities.

6.5.6 Conclusion

This validation study demonstrates that the ABI PRISM® 7000 Real-Time PCR system with SDS Software v1.2.3 is a robust, reliable, and reproducible system for performing DNA quantification using the Quantifiler kits.

When comparing 7000 System (SDS Software v1.2.3) results (C_T , slope, and R^2 values) to results obtained using the previously validated 7000 System (SDS Software v1.0):

- Small percentage differences are observed in C_T , slope, and R^2 values.
- Differences in calculated quantities are minimal for unknown samples using the 7000 System (SDS Software v1.2.3) and 7000 System (SDS Software v1.0).
- The differences observed should have little effect on resulting STR amplification based on calculated DNA quantities.
- No significant difference is observed between C_T values and calculated quantities derived using Auto Baseline and Manual analysis modes.

Bibliography

Afonina, I., Zivarts, M., Kutyavin, I., *et al.* 1997. Efficient priming of PCR with short oligonucleotides conjugated to a minor groove binder. *Nucleic Acids Res.* 25:2657–2660.

DNA Advisory Board. Quality assurance standards for forensic DNA testing laboratories (approved October 1998). Forensic Science Communications 2000;2(3). Available at:
<http://www.fbi.gov/programs/lab/fsc/backissu/july2000/codispre.htm>.

Förster, V.T. 1948. Zwischenmolekulare Energiewanderung und Fluoreszenz. *Ann. of Phys. (Leipzig)* 2:55–75.

Green, R., Roinestad, I., Boland, C., Hennessy, L., *Developmental Validation of the Quantifiler Real-Time PCR Kits for the Quantification of Human Nuclear DNA Samples*. Journal of Forensic Science, July, 2005, Vol 50, No. 4.

Kutyavin, I.V., Lukhtanov, E.A., Gamper, H.B., and Meyer, R.B. 1997. Oligonucleotides with conjugated dihydropyrroloindole tripeptides: base composition and backbone effects on hybridization. *Nucleic Acids Res.* 25:3718–3723.

Lakowicz, J.R. 1983. Energy Transfer. In *Principles of Fluorescence Spectroscopy*, New York: Plenum Press 303–339.

Martens, H. and Naes, T. 1989. *Multivariate Calibration*, Chichester: John Wiley & Sons.

Scientific Working Group on DNA Analysis Methods. Revised validation guidelines (approved July 2003). Available at:
http://www.fbi.gov/hq/lab/fsc/backissu/july2004/standards/2004_03_standards02.htm.

Index

Numerics

- 5' nuclease assay [1-3](#), [1-4](#)
- 7000 SDS
 - analysis settings, checking [5-2](#)
 - baseline settings [4-3](#)
 - detectors, creating [2-11](#)
 - detectors, settings [5-2](#)
 - fluorescence detection on [1-7](#)
 - PCR reactions, running on [3-7](#)
 - plate document, analyzing [4-3](#)
 - plate document, setting up [2-10](#)
 - results, viewing [4-4](#)
 - starting [2-5](#)
 - supported configuration [2-4](#)
 - template, setting up [2-22](#)
 - threshold settings [4-3](#)
 - validation SDS, v1.0 software [6-3](#)
 - validation, SDS v1.2.3 software [6-66](#)
- 7500 Real-Time PCR System, validation, SDS v1.2.3 software [6-51](#)
- 7900HT SDS
 - analysis settings, checking [5-3](#)
 - baseline settings [4-7](#)
 - fluorescence detection on [1-7](#)
 - PCR reactions, running on [3-9](#)
 - plate document, analyzing [4-7](#)
 - plate document, setting up [2-31](#)
 - results, viewing [4-8](#)
 - starting [2-27](#)
 - supported configuration [2-26](#)
 - template, setting up [2-40](#)
 - threshold settings [4-7](#)
 - validation, SDS v2.0 software [6-37](#)
- 9600 Emulation box
 - selection on the 7000 SDS [2-20](#)
 - selection on the 7900HT SDS [2-39](#)

A

- amplification plot
 - abnormal plots in one column, example of [5-14](#)
 - about [1-10](#)
 - baseline spikes, example of [5-13](#)
 - example [1-10](#)
 - inconsistent replicates, example of [5-12](#)
 - jagged plot, example of [5-13](#)
 - phases of [1-11](#)
 - print setup [4-6](#)
 - printing on the 7900HT SDS [4-10](#)
 - troubleshooting [5-12](#)
 - undefined plots, example of [5-14](#)
 - viewing on the 7000 SDS [4-5](#)
 - viewing on the 7900HT SDS [4-9](#)
- analysis settings
 - checking on the 7000 SDS [5-2](#)
 - checking on the 7900HT SDS [5-3](#)
 - verifying on the 7000 SDS [4-3](#)
 - verifying on the 7900HT SDS [4-7](#)
- Applied Biosystems
 - contacting [xiv](#)
 - customer feedback on documentation [xiv](#)
 - Information Development department [xiv](#)
 - Services and Support [xiv](#)
 - Technical Support [xiv](#)

B

- background, assay [6-35](#)
- baseline
 - about [1-13](#)
 - settings for the 7000 SDS [4-3](#)
 - settings for the 7900HT SDS [4-7](#)
- biohazard warning [xii](#)
- biohazardous waste, handling [xii](#)
- biological hazard safety. *See* biohazard warning
- bold text, when to use [vii](#)

C

- casework sample analysis 6-46
- CAUTION, description ix
- chemical safety guidelines x
- chemical waste
 - hazards xi
 - safety guidelines xi
- cleavage, in 5' nuclease assay 1-5
- clipped results, exporting (7900HT SDS) 4-10
- comparisons with other methods
 - on the 7000 SDS 6-27
 - on the 7900HT SDS 6-42
- composite spectrum 1-8
- computer, starting
 - for the 7000 SDS 2-5
 - for the 7900HT SDS 2-27
- contents of kit 1-16
- conventions
 - bold text vii
 - for describing menu commands vii
 - IMPORTANT! viii
 - in this guide vii
 - italic text vii
 - Notes viii
 - user attention words viii
- C_T. *See* threshold cycle
- customer feedback, on Applied Biosystems documents xiv

D

- DANGER, description ix
- data collection 1-8
- degraded DNA studies 6-23
- Detector Manager 2-14
- detectors (7000 SDS)
 - adding to plate document 2-14
 - analysis settings for 5-2
 - applying for standards 2-15
 - applying for unknown samples 2-17
 - creating 2-11

- detectors (7900HT SDS)
 - applying for standards 2-35
 - applying for unknown samples 2-36
 - copying to plate document 2-34
 - creating 2-32
- DNA quantification standards
 - dilution series, example 3-3
 - dilution series, guidelines for 3-2
 - guidelines for preparing 3-4
 - materials required to prepare 3-2
 - omitting Std. 8 5-5
 - preparing 3-4
 - reaction recommendation 3-6
 - See also* standards
- documentation, related xiv

E

- equipment, not included with Quantifiler kits 1-17
- exponential phase. *See* geometric phase
- exporting results
 - on the 7000 SDS 4-6
 - on the 7900HT SDS 4-10

F

- fluorescence, detection of 1-7

G

- geometric phase, amplification plot 1-11
- glycogen, adding to T₁₀E_{0.1} buffer 3-2
- guidelines
 - chemical safety x
 - chemical waste safety xi
 - waste disposal xii

H

- hazards
 - biological xii
 - chemical waste xi

I

- IMPORTANT, description [ix](#)
- Information Development department,
 - contacting [xiv](#)
- instrument, powering on
 - for the 7000 SDS [2-6](#)
 - for the 7900HT SDS [2-27](#)
- Internal PCR Control system, *See* IPC system
- IPC system
 - about assay [1-3](#)
 - components [5-10](#)
 - detectors for, creating on the 7000 SDS [2-13](#)
 - detectors for, creating on the 7900HT SDS [2-34](#)
 - interpreting results of [5-10](#)
 - invalid results from [5-11](#)
- italic text, when to use [vii](#)

L

- linear phase, amplification plot [1-11](#)

M

- manuals. *See* documentation, related
- materials, not included with Quantifiler kits [1-17](#)
- menu commands, conventions for
 - describing [vii](#)
- minor groove binder, description [1-4](#)
- mixture studies
 - on the 7000 SDS [6-21](#)
 - on the 7900HT SDS [6-40](#)
- MSDSs
 - description [x](#)
 - obtaining [x](#), [xiv](#)
 - referring to [x](#), [xi](#)
- multicomponent data
 - exporting (7900HT SDS) [4-10](#)
 - printing (7900HT SDS) [4-10](#)
 - processing of [1-9](#)

N

- negative results [5-11](#)
- New Document dialog box
 - 7000 SDS [2-10](#)
 - 7900HT SDS [2-31](#)
- nonfluorescent quencher, description [1-4](#)

O

- Optical Adhesive Cover, sealing plate
 - with [3-6](#)

P

- passive reference
 - multicomponent analysis, use in [1-9](#)
 - normalization using [1-10](#)
 - selecting in Well Inspector (7000 SDS) [2-15](#), [2-17](#)
 - selecting in Well Inspector (7900HT SDS) [2-36](#), [2-37](#)
- PCR
 - inhibition [5-11](#)
 - kinetic analysis of [1-10](#)
 - phases of [1-11](#)
 - process in 5' nuclease assay [1-4](#)
 - reactions, preparing [3-5](#)
 - reactions, running on the 7000 SDS [3-7](#)
 - reactions, running on the 7900HT SDS [3-9](#)
 - standard, preparing [3-2](#)
- plate document (7000 SDS)
 - analyzing [4-3](#)
 - creating blank [2-10](#)
 - creating from a template [2-24](#)
 - detectors, adding [2-14](#)
 - detectors, applying for standards [2-15](#)
 - detectors, applying for unknown samples [2-17](#)
 - detectors, creating [2-11](#)
 - how used [2-7](#)
 - sample names, adding [2-18](#)
 - saving [2-21](#)
 - setup examples [2-8](#)
 - template, creating [2-22](#)
 - template, setting up [2-22](#)
 - thermal cycler conditions, setting [2-19](#)
 - types [2-7](#)

plate document (7900HT SDS)
 analyzing 4-7
 creating blank 2-31
 creating from a template 2-42
 detectors, applying for standards 2-35
 detectors, applying for unknown
 samples 2-36
 detectors, copying 2-34
 detectors, creating 2-32
 how used 2-28
 sample names, adding 2-37
 saving 2-39
 setup examples 2-29
 template, creating 2-40
 template, setting up 2-40
 thermal cycler conditions, setting 2-38
 types 2-28

plateau phase, amplification plot 1-11

polymerization, in 5' nuclease assay
 completion of 1-6
 process 1-5

precision
 on the 7000 SDS 6-4
 on the 7900HT SDS 6-37

printing results
 on the 7000 SDS 4-6
 on the 7900HT SDS 4-10

probes, about 1-4

pure dye spectra, how used 1-9

Q

Quantifiler Human detector
 creating (7000 SDS) 2-12
 creating (7900HT SDS) 2-32

Quantifiler Y detector
 creating (7000 SDS) 2-12
 creating (7900HT SDS) 2-33

quantity
 assessing 5-16
 entering for detector (7000 SDS) 2-16
 exporting (7000 SDS) 4-6
 printing (7000 SDS) 4-6
 viewing (7000 SDS) 4-5
 viewing (7900HT SDS) 4-9

R

R² value
 interpreting 5-4
 viewing (7000 SDS) 4-4
 viewing (7900HT SDS) 4-8

radioactive waste, handling xii

raw data
 about 1-8
 printing (7000 SDS) 4-6
 printing (7900HT SDS) 4-10

reactions
 examples of arranging 2-8, 2-29

real-time data analysis 1-10

regression line formula 5-4

replicates
 grouping in exported results (7900HT
 SDS) 4-10
 results, how grouped 2-18, 2-37

report
 exporting and printing on the 7000
 SDS 4-6
 exporting and printing on the 7900HT
 SDS 4-10
 viewing on the 7000 SDS 4-5

reporter signal, normalized
 about 1-10
 viewing in amplification plot (7000
 SDS) 4-4
 viewing in amplification plot (7900HT
 SDS) 4-8

reproducibility 6-7

results
 exporting (7900HT SDS) 4-10
 viewing (7000 SDS) 4-4
 viewing (7900HT SDS) 4-8

R_n. See reporter signal, normalized

ROX
 selecting in Well Inspector (7000
 SDS) 2-15, 2-17
 selecting in Well Inspector (7900HT
 SDS) 2-36, 2-37

S

safety

- biological hazards
- chemical waste [xi](#)

safety alert words

- CAUTIONS [ix](#)
- DANGERS [ix](#)
- description [ix](#)
- IMPORTANTS [ix](#)
- WARNINGS [ix](#)

sample name

- adding for unknown samples (7000 SDS) [2-18](#)
- adding for unknown samples (7900HT SDS) [2-37](#)
- entering for detector (7000 SDS) [2-16](#)
- entering for detector (7900HT SDS) [2-36](#)

SDS document. *See* *.sds file

*.sds file

- description (7000 SDS) [2-7](#)
- description (7900HT SDS) [2-28](#)
- saving (7000 SDS) [2-21](#)
- saving (7900HT SDS) [2-39](#)
- See also* plate document (7000 SDS) and plate document (7900HT SDS)

SDS software

- starting for 7000 SDS [2-6](#)
- starting for 7900HT SDS [2-27](#)

SDS template. *See* *.sdt file

*.sdt file

- description (7000 SDS) [2-7](#)
- description (7900HT SDS) [2-28](#)
- saving (7000 SDS) [2-23](#)
- See also* template

sensitivity

- of assay [5-16](#)
- tests [6-16](#)

Services and Support, obtaining [xiv](#)

setup table, exporting (7900HT SDS) [4-10](#)

single plate document. *See* *.sds file (7900HT SDS)

slope of standard curve

- about [5-4](#)
- differs significantly from -3.33 [5-6](#)
- interpreting [5-5](#)

viewing (7000 SDS) [4-4](#)

viewing (7900HT SDS) [4-8](#)

specificity

- with bacterial pools panel [6-14](#)
- with human DNA panel [6-10](#)
- with non-human panel [6-11](#)

stability [6-17](#)

standard curve

- about results [5-4](#)
- differences in C_T values of replicates [5-9](#)
- interpreting [5-4](#)
- outlier in, example of [5-8](#)
- print setup (7000 SDS) [4-6](#)
- replicates, example of four [5-9](#)
- straight horizontal line [5-7](#)
- troubleshooting [5-6](#)
- viewing (7000 SDS) [4-4](#)
- viewing (7900HT SDS) [4-8](#)

standards

- applying detectors for (7000 SDS) [2-15](#)
- applying detectors for (7900HT SDS) [2-35](#)
- See also* DNA quantification standards

stochastic effects [5-16](#)

storage recommendations, for kits [1-16](#)

strand displacement, in 5' nuclease assay [1-5](#)

T

T₁₀E_{0.1} buffer, contents of [3-2](#)

TaqMan[®] MGB probe [1-4](#)

targets, about [1-3](#)

task

- selecting for detector (7000 SDS) [2-16](#)
- selecting for detector (7900HT SDS) [2-35](#)

Technical Support, contacting [xiv](#)

template

- creating plate document from (7000 SDS) [2-24](#)
- creating plate document from (7900HT SDS) [2-42](#)
- saving (7000 SDS) [2-23](#)
- saving (7900HT SDS) [2-41](#)
- setting up for plate documents (7000 SDS) [2-22](#)
- setting up for plate documents (7900HT SDS) [2-40](#)

template plate document. See *.sdt file

text conventions [vii](#)

thermal cycler conditions

- printing (7900HT SDS) [4-10](#)
- setting (7000 SDS) [2-19](#)
- setting (7900HT SDS) [2-38](#)

threshold

- settings for the 7000 SDS [4-3](#)
- settings for the 7900HT SDS [4-7](#)

threshold cycle

- calculation of [1-13](#)
- exporting (7000 SDS) [4-6](#)
- in standard curve [5-4](#)
- normal range for IPC system [5-11](#)
- printing (7000 SDS) [4-6](#)
- relationship to initial template amount [1-14](#)
- viewing in amplification plot (7000 SDS) [4-4](#)
- viewing in amplification plot (7900HT SDS) [4-8](#)
- viewing in standard curve (7000 SDS) [4-4](#)
- viewing in standard curve (7900HT SDS) [4-8](#)

training, information on [xiv](#)

U

unknown samples

- detectors, applying for (7000 SDS) [2-17](#)
- detectors, applying for (7900HT SDS) [2-36](#)
- names, adding (7000 SDS) [2-18](#)
- names, adding (7900HT SDS) [2-37](#)

user attention words, described [viii](#)

V

validation

- 7000 SDS, SDS v1.0 software [6-3](#)
- 7000 SDS, SDS v1.2.3 software [6-66](#)
- 7500 Real-Time PCR System, SDS v1.2.3 software [6-51](#)
- 7900HT Real-Time PCR System, SDS v2.0 software [6-37](#)
- background, assay [6-35](#)
- comparisons with other methods [6-27](#)
- comparisons with other methods (7900HT SDS) [6-42](#)
- degraded DNA studies [6-23](#)
- importance of [6-2](#)
- mixture studies [6-21](#)
- mixture studies (7900HT SDS) [6-40](#)
- precision [6-4](#)
- precision (7900HT SDS) [6-37](#)
- reproducibility [6-7](#)
- sensitivity [6-16](#)
- specificity with bacterial pools panel [6-14](#)
- specificity with human DNA panel [6-10](#)
- specificity with non-human panel [6-11](#)
- stability [6-17](#)

W

WARNING, description [ix](#)

waste disposal, guidelines [xii](#)

Y

y-intercept of standard curve

- interpreting [5-4](#)
- viewing (7000 SDS) [4-4](#)
- viewing (7900HT SDS) [4-8](#)

Headquarters

850 Lincoln Centre Drive
Foster City, CA 94404 USA
Phone: +1 650.638.5800
Toll Free (In North America): +1 800.345.5224
Fax: +1 650.638.5884

Worldwide Sales and Support

Applied Biosystems vast distribution and service network, composed of highly trained support and applications personnel, reaches 150 countries on six continents. For sales office locations and technical support, please call our local office or refer to our Web site at www.appliedbiosystems.com.

www.appliedbiosystems.com



Applera Corporation is committed to providing the world's leading technology and information for life scientists. Applera Corporation consists of the Applied Biosystems and Celera Genomics businesses.

Printed in the USA, 04/2006
Part Number 4344790 Rev. D

an **Applera** business