

Hopkins® 2-in-1 Digital Baby Scale User's Manual


#688330

Technical Specifications

- Operating Temperature Range: 50°F ~ 104°F (10° ~ 40°C)
- Environmental Humidity: < 80%
- Maximum Weight Capacity: 44lb/20kg
- Graduation: 0.01lb/5g
0-22lbs (0-10kg) ± 0.5oz (15g)
22-44lbs (10-20kg) ± 1oz (30g)
- Power Supply: 4 AA Batteries (INCLUDED), Automatic Shut Off
- Background Color: Blue Backlight
- Low Power Warning Display: **LO**
- Error Warning Display: **ERR**
- Over Load Warning Display: **OL**
- Units: **KG** or **LB**
- Memory Function: Save and Recall Previous Weight Measurement
- Hopkins® 3 Year Warranty Guarantee - please register your scale at: www.hmponline.com/registration

Operation

- After the batteries are inserted, press the **“ON/TARE”** key to start the scale until it self-zeroes. **“0.000”** will appear when the scale is ready for weighing.
- Gently place the baby on the tray, the weight will be shown in the window after a few seconds. If the baby remains on the tray after its weight value has shown, the weight value will continue to flash for 5 minutes and then turn off automatically.
- If the baby is removed from the tray immediately after the weight has been displayed, the scale will return to **“0.000”** and then turn off automatically within 5 minutes.
- To save the current weight into memory, press and hold the **“MEM”** key when the weight value is stable. The scale will display **“SET”** when the value has been saved.
- To recall the prior saved weight, press the **“MEM”** key while the scale is zeroed **“0.000”**.
- To weigh a toddler or child under 44 lbs, remove the baby weighing tray by gently squeezing the two release handles underneath the scale to separate the tray from the scale. To reattach the tray, simply align the tray and press down gently to snap into place.

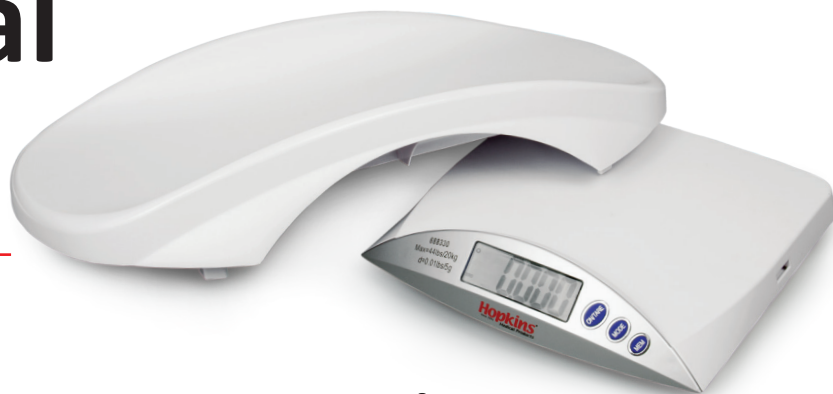

Button Functions

- **ON/TARE Button:** Press to power on the scale and start up weighing process. To zero or calibrate the scale, press the key again to clear the display and **“ZERO”** will show in the upper left-hand corner of the LCD. To zero with an object press **“ON/TARE”** with object (i.e. blanket or scale liners) Hold this key for 3 seconds to turn off the scale.
- **MODE Button:** Press the key to toggle between kg and lb. The unit of measure indicator is shown on the right side of the LCD as kg or lb.
- **MEM Button:** While the scale displays **“0.000”**, press the **“MEM”** to recall the previous saved weighing value. A stable value is indicated by the **“0”** in the top left corner of the display. Pressing the **“MEM”** key during a stable value will save the weight, an unstable value will not be saved.

Precautions

- The scale should be placed on a hard flat surface for accurate measurement and patient safety.
- Replace batteries immediately when the screen displays **“LO”**.
- Keep the scale clean to ensure proper functionality.
- Remove the batteries for long term storage. Battery corrosion is not covered under warranty.
- Place the patient in the center of the tray to ensure weighing accuracy.
- Overloading the scale with weight heavier than 45lb / 23kg is not recommended and can damage the scale, voiding the warranty.

Hopkins®
Since 1945
Medical Products