

no!
no!
HAIR +

User Manual



Contents

Safety Instructions	4
Meet no!no! Hair Plus.....	6
Before you no!no! with no!no! Hair Plus	8
Working with no!no! Hair Plus	10
Thermicon™ Tips.....	16
Treating Facial Hair	18
Maintenance	19
Troubleshooting	20
Technical Information	21

Safety Instructions

READ ALL INSTRUCTIONS BEFORE USE.

ALL WARNINGS AND CAUTIONS SHOULD BE STRICTLY FOLLOWED

As with any electrical device, certain precautions must be taken in order to ensure your safety.

WARNINGS – To reduce the risk of burns or injury:

- Do not use this device on the genitals, nipple area, upper lip or around the eye.
- On concave, folded or loose skinned areas, stretch the skin tight and make sure it is as straight as possible to avoid cuts and burns.
- Do not use this device on sunburned skin, open cuts, bruises, blisters, moles or on raised areas such as scars.
- If you have any skin condition other than those listed above, like skin cancer or cold sores, please consult your physician before using this device.
- Before using no!no! Hair Plus make sure your skin surface is completely dry.
- Keep out of reach of children.
- Do not use deodorants, lubricants, liquids, creams, gels, water, and etc. before use.
- Do not allow the thermodynamic wire to touch your skin.
- Do not use the Thermicon™ Tip if the wire is bent or damaged in any way.
- Only use this device for its intended purpose as described in this manual.

CAUTIONS - to reduce the risk of fire, electric shock, etc.

- Do not reach for a device that has fallen into water while it is plugged in. Unplug the power supply before retrieving it from the water.
- Do not place or store device near a sink or bathtub.




- This device should never be left unattended when plugged in.
- Turn off the device while changing the Thermicon™ Tip.
- Do not use the power supply if the cord or plug is damaged in any way.
- Do not use this device near or rest it on flammable objects such as paper or cloth.
- Always turn off and unplug this device when you have finished working with it.
- Always turn off and unplug the power cord from the electrical outlet before cleaning or maintenance.
- Do not place this device near water or other liquids.
- This device should only be used with the supplied power supply.
- Do not open the outer casing of the unit. This might harm the unit and it will void any warranty.

Meet no!no! Hair Plus

Forget everything you've ever learned about hair removal. no!no! Hair Plus is about to show you something new and exciting – Painless hair removal with professional salon-like results, all in the comfort and convenience of your home.

no!no! Hair Plus is a professional treatment that, when used continuously, may lead to a reduction in hair regrowth. This is science, not magic. Results take time, patience and a strong commitment to the treatment's procedures. And since every body is unique, results will vary. The average user should notice results in 6-8 weeks, for others it might take a bit longer.

For the safest, most effective treatment, read the entire user manual before you begin to no!no! with your no!no! Hair Plus.

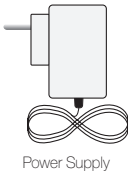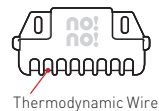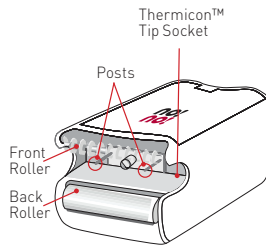

Before you no!no! with no!no! Hair Plus

Get the best results and performance by learning as much as you can before you no!no! with your no!no! Hair Plus. Below is a brief description about what to expect and what to look out for.

- **Test run** – Before you use no!no! Hair Plus for the first time, and before treating new areas, always test no!no! Hair Plus on a small section of the area you want to treat using the lowest treatment level.
- **Track your progress!** Before you use no!no! Hair Plus for the first time, study the area you want to treat. Then, once every 2 weeks, check your progress. But remember, it takes 6-8 weeks for initial results to be noticed*.
- **The right length** – Hair should be about stubble length, between 1-3 mm, for best treatment.
- **Aroma therapy** – no!no! Hair Plus burns hair and burnt hair smells. But that's how you know its working. Try aromatic candles or other air fresheners. An after treatment cream or moisturizer like no!no! Smooth will help mask the smell even further.
- **That feels weird** - Some people experience a mild tingling sensation when they first use no!no! Hair Plus, especially around sensitive areas. This should fade then disappear after the first few treatments.
- **For the sensitive type** – If your skin feels uncomfortable during treatment, reduce the treatment level and/or take short breaks after repeatedly gliding on any one area. If it continues, stop treatment completely.
- **Keep it flat** – On concave, folded or loose skinned areas, stretch the skin and make sure it is as straight as possible to avoid cuts and burns.


- **Remove stubble** – Buff after you use no!no! Hair Plus for a smooth, polished look and feel.
- **Keep it smooth** – When you finish your treatment, apply an alcohol-free lotion or cream, like no!no! Smooth, to keep your skin moisturized.
- **Patience is a virtue** - Initial results may take 6-8 weeks* to appear.
- **The first 6-8 weeks** - For best results, use no!no! Hair Plus 2 -3 times a week for the first 6-8 weeks. You don't have to wait for a lot of hair to grow back. Just 1mm of stubble is enough to get the smooth results you are looking for.
- **After the first 6-8 weeks** – You can reduce the treatments to once or twice a week and as time goes on, just use no!no! Hair Plus as needed.
- **Check the Thermicon™ Tips** - Always make sure the thermodynamic wire is completely intact before use. If it is bent or torn, replace it before your next treatment.
- **Clean the Thermicon™ Tips** – Use the cleaning brush frequently to remove hair debris, especially when working on long hair. This will help prolong the life of the Thermicon™ Tips.

*Results may vary.

Working with no!no! Hair Plus


Get Ready

1. Remove the safety seal.

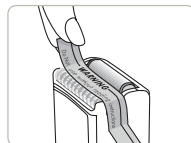

2. Insert the Thermicon™ Tip


- Pick up the device with the rollers facing you. Inside the Thermicon™ Tip Socket you will see 2 posts.
- Hold the Thermicon™ Tip so you can see the 2 holes

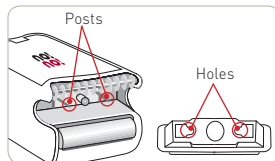

3. Line up the posts in the socket with the holes in the Thermicon™ Tip. Make sure the no!no! logos are facing the same direction. Gently push into place. Do not use force.



4. Be sure the Thermicon™ Tip is correctly in place.

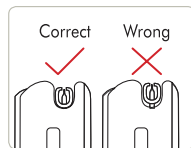

NOTE:

To remove the Thermicon™ Tip

- Turn the device off.
- Grasp it by both ends.
- Gently pull.

Get Set

Prepare your skin

no!no! Hair Plus must be used on clean, dry skin. Always be sure to clean your skin of any gels, oils, lotions or deodorants before you begin to use it.


■ Plug it in

Connect the power supply to your no!no! Hair Plus and plug it into an electrical outlet.




■ Turn it on

Press the on/off button found on the no!no! Hair Plus power supply.

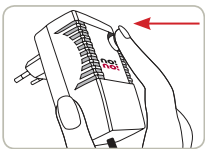


■ Adjust the Treatment Level

For your comfort, there are 3 different treatment levels.

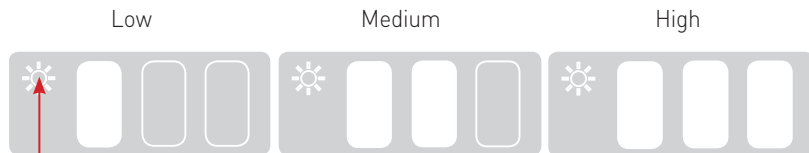

Treatment Level Icon

Select the level by pressing the button on top of your no!no! Hair Plus until the desired treatment level appears.

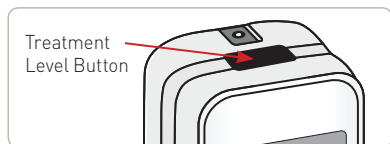


- Press once for Medium, twice for High and once more to return to Low.
 - When first using no!no! Hair Plus or when starting to treat a new area, test it on a small section of skin using the lowest treatment level.
 - If you can, raise the level and work at the highest level possible while still maintaining comfort.
 - If your skin feels hot or uncomfortable at anytime, stop immediately. Lower the treatment level and try again.
- * When you turn on no!no! Hair Plus again, it will automatically be set to the last programmed treatment level.

Master the Slow Glide

If you are treating facial hair, please see the section - Treating Facial Hair.

Crucial to the treatment is the Slow Glide. To master this technique:

1. Hold no!no! Hair Plus at a 90° angle with the rollers flat against the skin and the status screen facing you.

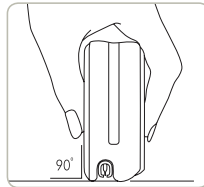

2. Gently glide no!no! Hair Plus over each section you want to treat using a slow, smooth, even motion.
 - Always glide in one direction, with the status screen facing you.
 - Do not stop in one place.
 - Glide over the section you are treating a few times until you feel it is smooth or nearly smooth. Don't over do it.



Master the Buffer

Buffing is the essential final step to a no!no! Hair Plus treatment. Buffing exfoliates your skin, removes crystallized hair and leaves your skin feeling silky smooth To Buff:

1. Hold the Buffer flaps
2. Firmly rub the buffer in either a circular or back and forth motion over the treated area.



Important Note:

Do not be alarmed when the number of bars changes on the Status Screen.


When you first start to Glide, the Status Screen will automatically switch to the Treatment Level display (1) to the Thermicon™ Tip display (2) when you begin to glide.



Check the icon to determine what information is being displayed.

Rounded areas

To work on more difficult areas, go slowly and follow the natural curve of the body. We recommend mastering the Glide on arms and legs before moving on to the knees, elbows, face and bikini line.

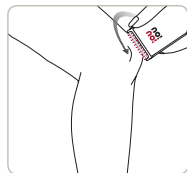

NOTE: The face and bikini line may be extra-sensitive to temperature. If it becomes uncomfortably warm during treatment, lower the treatment level or work on a different area for a while. **Do not proceed to the genital area.**

Watch the Light


The signal light by the Thermicon™ Tip will guide you throughout your no!no! Hair Plus experience. It lets you know when you are working at the right speed and warns you if the Thermicon™ Tip is torn, damaged or no longer working. The Signal Light is visible only during the Slow Glide part of treatment.

When working at the right speed, the light will remain steady. If the light flickers, adjust your speed accordingly. If the light on top of the no!no! Hair Plus turns red – check the Thermicon™ Tip.

Blue	Red	What it means
On and steady	Off	Correct gliding speed
Off	On	No Thermicon™ Tip in socket
Off	On	Depleted Thermicon™ Tip
Off	On	Damaged Thermicon™ Tip
Flickering	Off	Incorrect gliding
Off	Off	No movement

Thermicon™ Tips

To avoid possible malfunction or harm, make sure to gently push the Thermicon™ Tip all the way into the socket.

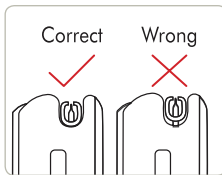


WARNING! Always check the wire in the Thermicon™ Tip before use. Never use a Thermicon™ Tip if the wire is bent, torn or damaged in any way.



Bent thermodynamic wire

© 2011 ThermoFisher Scientific

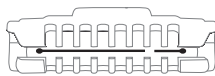





Torn thermodynamic wire

The Thermicon™ Tip Status Screen:

The Thermicon™ Tip Status Screen automatically appears as you begin to Glide and is indicated by the appearance of the Tip icon.

Replace the Thermicon™ Tip...

- When there are no bars on the status screen
- When the Thermicon™ Tip icon flickers
- If it is no longer working efficiently
- If the wire is damaged in any way

Thermicon™ Tip Status	
New	Slightly Used
 Thermicon™ Tip Icon	
Used	Replace
	

Treating Facial Hair

Your face is very sensitive. Take a few minutes and read these important tips for the safest, most effective treatment.

And remember, **If there is any discomfort – STOP.**

■ Test

Before you begin, use the instructions below and treat a small, discreet section of your face to see how your skin reacts to a no!no! Hair Plus facial hair treatment. Remember to always start with the lowest treatment level and work your way up.

If there is any discomfort – STOP.

■ Chin & Jaw Line

Working from just below the jaw or chin line, glide your no!no! up and around the curve.

If you need to, pull skin straight and/or purse your lips inward to give your chin and jaw a straighter, smoother line.

If there is any discomfort – STOP.

■ Sideburns

Start at the base of the chin and glide upwards following the curve of your face. If you need to, pull your skin tight and straight for a flat, smooth working surface.

If there is any discomfort – STOP.

Important Note:

The areas just below the eye, around the eyelid and the upper lip are especially sensitive.

DO NOT treat these areas.

Maintenance


Do not use any liquids of any kind to clean the Thermicon™ Tip or rollers.

Thermicon™ Tips

It is important to keep the Thermicon™ Tip free of hair debris. Allowing hair to build up may cause it to burn out prematurely or malfunction. Check and clean the tip periodically during each treatment session.

To clean the Thermicon™ Tip:

1. Turn off your no!no! Hair Plus.
2. Remove the Thermicon™ Tip.
3. Gently brush clean.

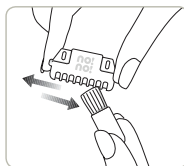


Rollers

Although not necessary to the performance of no!no! Hair Plus, keeping the rollers clean is all part of proper maintenance.

Use the cleaning brush to remove excess hair and debris from the rollers.

Brush debris out and away from the socket in order to prevent any particles from falling into the socket.

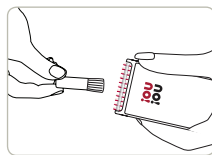

If any hair or debris fall into the socket, this may cause no!no! to not function properly. Gently blow any particles out of the socket before beginning the next treatment.

Troubleshooting

If no!no! Hair Plus does not start when you glide it on the skin:

- Press the On/Off button on the power supply again to make sure it is turned on.
- Check the Status Screen while gliding. If no bars appear and the Red signal light is on, replace the Thermicon™ Tip.
- Make sure the Thermicon™ Tip is not bent or torn. Replace it if necessary.

If the device is not working properly, contact your point of purchase for examination or repair. The no!no! Hair Plus warranty is only valid if no!no! Hair Plus has been used according to the instructions in this guide. Retain all receipts and original packaging.

Technical Information

Technical Specifications:

Weight: 130gr

Physical Dimensions: 3.8 X 3.3 X 81mm
(1.5 X 1.3 X 3.2 in)

Mode of Operation: Continuous

Operating Conditions:

Temperature: 5° C-30° C

Humidity: 30%-80% RH

Rating: 100-240VAC; 50-60 HZ; 0.3A

Transport and Storage

Temperature: -20° C to 80° C

Humidity: 05-95% RH

Atmospheric pressure: 0.1-1.5 atm.

Electrical Safety:

Power Supply:  Class II equipment – double insulation

This appliance is not suitable for use in the presence of flammable anesthetic mixture with air or with Oxygen or Nitrous Oxide.

The Quality Management System of Radiancy Ltd. complies with the Quality Management Standard ISO 13485.

no!no! devices comply with IEC 60335 and with Low Voltage Directive (Electric Safety) LVD: 2006/95/EC

Design and technical specifications are subject to change without notice.

Manufacturer:

Radiancy Israel LTD.

5 Hanagar Street

45240 Hod Hasharon

Israel

EC Representative:

Obelis S.A Av.de Tervuren 34 Bte. 44

B-1040 Brussels, BELGIUM

Tel: 32 (0) 2 732 59 54

Fax: 32 (0) 2 732 60 03

GSM: 07545 46 60

Email: mail@obelis.net





no!
no!
HAIR+