


# Conductivity Monitor CM-P



## User Manual

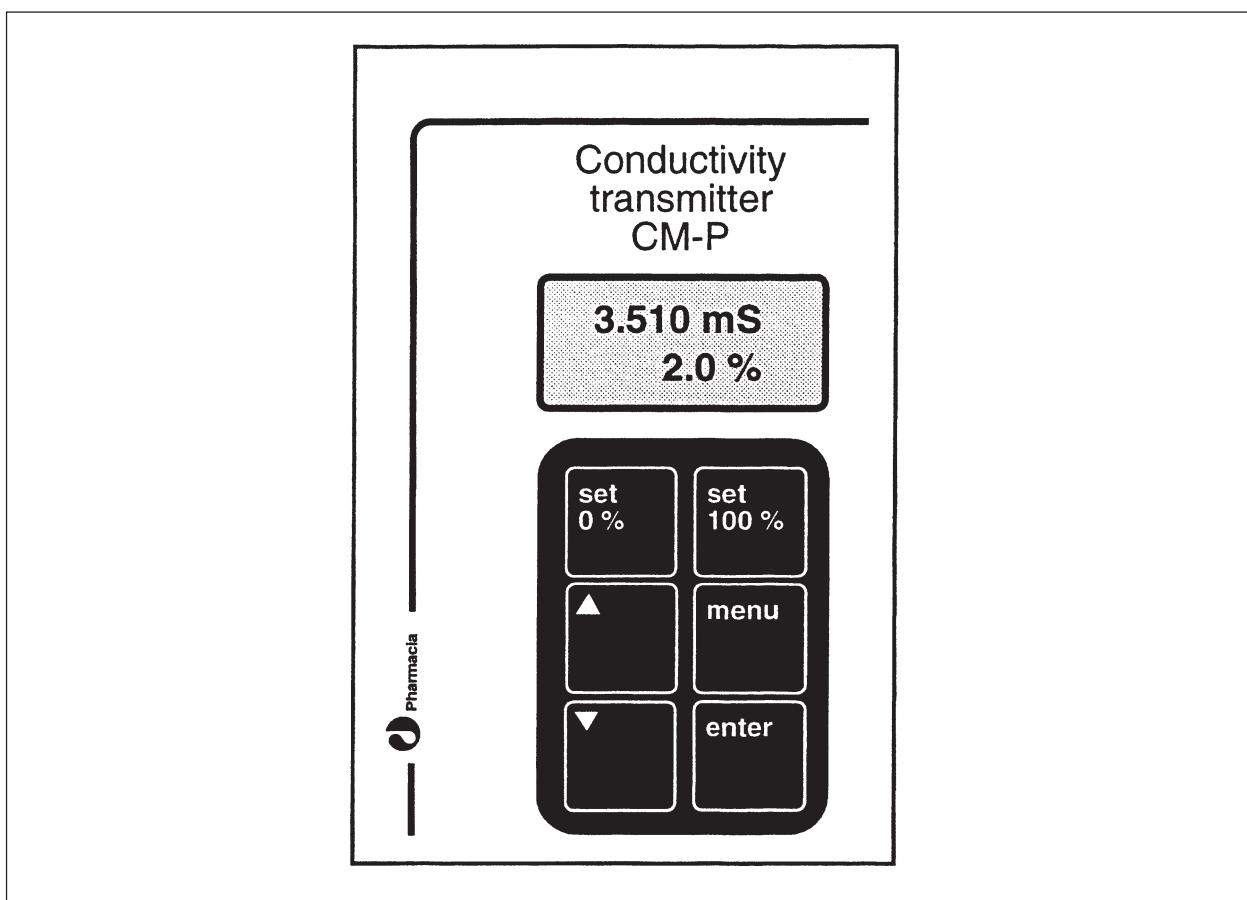

# Contents

1. Introduction .....	3
2. General Description .....	4
2.1 Principle of operation .....	9
2.2 Transmitter .....	5
2.2.1 Front panel .....	5
2.3 Flow cell .....	6
3. Technical Specifications .....	7
4. Installation .....	8
4.1 Connection of the flow cell .....	8
4.2 Transmitter set up .....	8
5. Operation .....	10
5.1 Setting the scale with buffers .....	10
5.2 Menu .....	11
5.3 Fixed scale setting .....	12
5.4 Temperature compensation .....	13
5.5 Summary of commands in menu .....	15
6. Maintenance .....	16
6.1 Cleaning the flow cell .....	16
6.2 Calibration mode .....	16
6.3 Calibration of the temperature sensor .....	17
6.4 Calibration of the cell constant .....	18
6.5 Replacing the flow cell .....	20
6.6 Summary of commands in Calibration Mode .....	21
7. Trouble Shooting .....	22
8. Accessories .....	23



# 1. Introduction

The Amersham Biosciences Conductivity Monitor, CM-P, is a high quality, high precision on-line monitor for conductivity measurements and gradient determination in liquid chromatography applications.

The monitor consist of two units, the flow cell and the transmitter. The flow cell contains the titanium electrodes used for conductivity measurement as well as a temperature sensor. The small size of the flow call allows it to be posiotioned anywhere in the liquid path, provided that the pressure specification for the cell is not exceeded.

The transmitter carries the six programming keys for the monitor. An easy to read 16 digit alphanumeric display shows the conductivity measured in Siemens/cm (mS/cm).

The transmitter communicates with the UNICORN controller over a serial interface using optomux protocol. Through this communication the conductivity measurement is transmitted to UNICORN. The electrical outputs are only used for a possible optional recorder.

The monitor has a very large dynamic range from 0.001 mS/cm to 999.9 mS/cm and therefore suitable for a wide range of applications, ranging from reversed phase to hydrophobic interaction chromatography. A unique design of the electronics is the microprocessor control of frequency variations within the working range. This gives excellent linearity and eliminates the need for range settings.

## 2. General description

### 2.1 Principle of operation

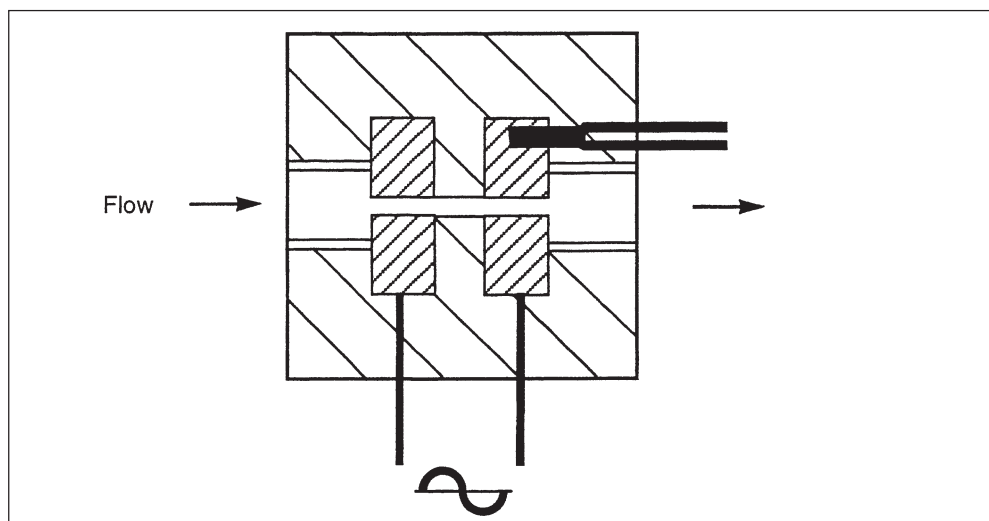


Fig. 1. The flow cell

The flow cell has two cylindrical shaped titanium electrodes positioned in the flow path of the cell. An alternating voltage is applied between the electrodes and the resulting current is measured and used to calculate the conductivity of the eluent. The microprocessor controls the AC frequency and automatically increases it with increasing conductivity between 50 Hz and 50 kHz. This unique design gives maximum linearity and true conductivity values.

The monitor measures conductivity over the complete working range. No range settings are required. However, to obtain a suitable output signal to the recorder and to get a relative conductivity value, it is possible to set a 0 to 100% range. This can be done with the buffers used or by selecting any fixed range between 0.001 mS/cm and 999.9 mS/cm.

The monitor displays the conductivity in mS. The conductivity is automatically calculated by multiplying the measured conductance by the flow cell's cell constant. The cell constant is pre-calibrated on delivery but can be measured with a separate calibration procedure.

One of the electrodes has a small temperature sensor for measuring the temperature of the eluent in the flow path. Temperature variations influence the conductivity and in some applications, when highly precise conductivity values are required, it is possible to program a temperature compensation factor that recalculates the conductivity to a set reference temperature.

## 2.2 Transmitter

The microprocessor-based transmitter communicates with the UNICORN controller over a serial interface using the optomux communication protocol. The transmitter fits into a 19" industry standard rack.

### 2.1.1 Front panel

The front panel consists of an alphanumeric LCD display window and a keyboard with six membrane keys.

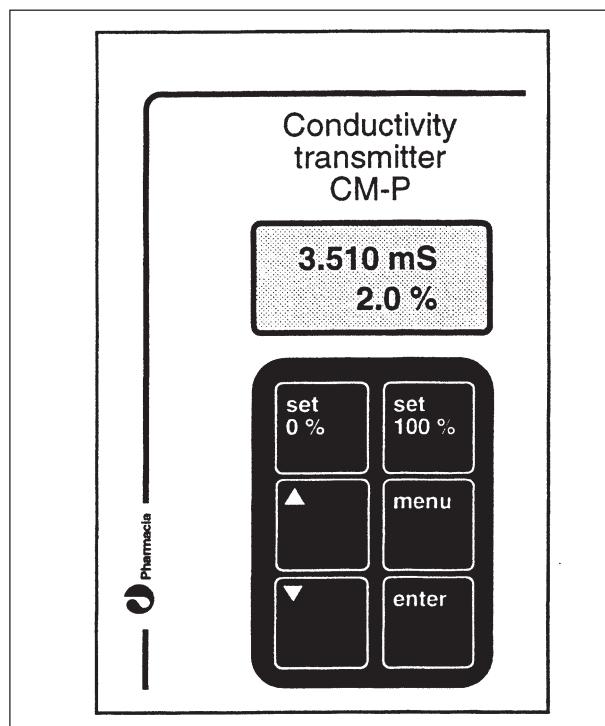


Fig. 2. The front panel

### Display

mS Displays the actual conductivity value in mS, (mS/cm) and  
% the relative conductivity in %.

### Keyboard



This key sets the 0% conductivity level with the low conductivity buffer. The set 0% key is pressed for 1 second while the low conductivity buffer is pumped through or injected into the flow cell. The voltage output to the recorder is set to 0 V for the measured conductivity.



This key sets the 100% conductivity level with the high conductivity buffer. The set 100% key is pressed for 1 second while the high conductivity buffer is pumped through or injected into the flow cell. The voltage output to the recorder is set to 1 V for the measured conductivity.



The menu is used to step between the functions below:

1. Set 100%. With this function any fixed conductivity value in mS/cm can be set for the 100% level.
2. Set 0%. With this function any fixed conductivity value in mS/cm can be set for the 0% value.

3. Set TC. A temperature compensation factor for the conductivity can be set. This can be necessary for high precision measurements. The temperature compensation factor is expressed in percentage increase of conductivity per °C increase in temperature. If the temperature compensation factor is unknown, a general approximate value for many common salt buffers used in chromatography applications.  
When the temperature compensation is set to 0, it is inactive.
4. Set TRcf. The temperature compensation factor must be realloted and adjusted to a set reference temperature in range 4-40 °C.



The "arrow down" and "arrow up" keys are used to select numerical values in the menu settings. The keys can be held down for a faster scroll.



The enter key is used to enter programmed values and to exit to Run Mode.

## 2.3 Flow cell

The flow cell exterior is made of PEEK and the interior between the titanium electrodes of fluoroplastics. Both ends of the flow cell have female M6 fittings. The signal cable which connects to the control unit has a 9 pole DSUB connector.



# 3. Technical Specifications

Operating mode	Auto range conductivity in mS/cm
Sensitivity range	0.001 mS/cm-999.9 mS/cm
Accuracy	2% full scale calibrated range or 0.01 mS/cm whichever is greater in the range 0.001-300 mS/cm
Noise	0.5% full scale calibrated range
Flow rate range	0-100 ml/min
Analogue output	0-1 V or 4-20 mA
Power supply	24 V DC
Operating conditions	4-40 °C, 20-95% RH, Atmospheric pressure 840-1060 mb, (84-106 kPa)
Power consumption	2 W
<i>Flow cell</i>	
Type	Flow through
Cell constant	50 cm <sup>-1</sup> , nominal
Internal volume	14 lll
Maximum operating pressure	5 MPa, 50 bar, 750 psi
Wetted parts	Fluoroplastics, titanium
Fittings	M-6 (6 mm metric)

## 4. Installation

### 4.1 Connection of the flow cell

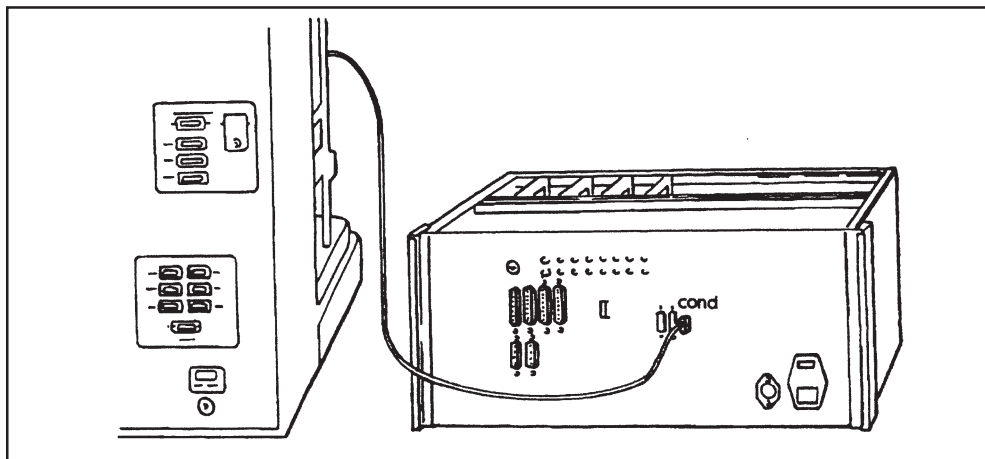


Fig. 3. Connection of flow cell to UNICORN control box.

1. Connect the signal cable from the cell to the socket "Cond" on the UNICORN Control box.

The transmitter is correctly set up on delivery of the system. Set-up parameters are saved in non-volatile memory, and do not need to be changed.

Make a note of the optomux communication parameters for your transmitter. If a faulty transmitter is replaced with a new unit, the optomux parameters for the new unit should be set to the same values as for the old transmitter.

To set parameters, start the monitor in start-up mode by pressing and simultaneously as the power is switched on, and keeping the buttons pressed for about 5 seconds. The display is dark for 5 seconds, then shows the program version for 5 seconds followed by the current optomux address, e.g.:

```
Optomux
addr 002
```

The address is a number between 0 and 255. Refer to your system documentation for the correct address. Press the arrow keys to change the address until the correct value is displayed, then press enter. The display now shows the communication rate, e.g.:

```
Baudrate
38400
```

Permissible baud rates are 1200, 2400, 4800, 9600, 19200 and 38400. Normally 38 400 is used. the correct setting is determined by the optopux

communication port setting in the controller. Use the arrow keys to choose the correct communication baud rate, then press enter. The display now shows the accumulated operation time and the number of starts, e.g.:

```
00021 h
```

```
0003 st
```

Press enter to continue.

# 5. Operation

## 5.1 Setting the with buffer

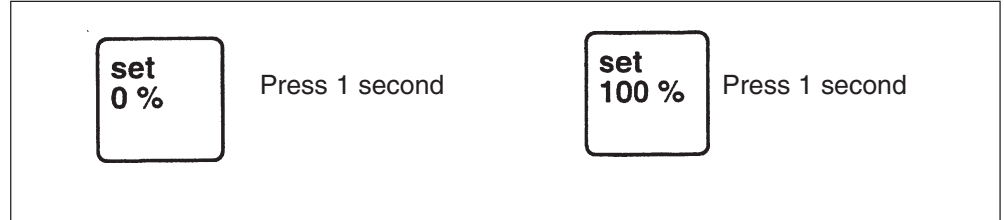




Fig. 4. Setting the scale.


When the conductivity monitor is switched on, it is in Run Mode and will display the measured conductivity. To set the output range to the recorder, the conductivity levels 100% and 0% can easily be set with the high and low conductivity buffers used for the gradient

1. Start the flow in the system with the high conductivity buffer.
2. When the conductivity reading stabilises, press  for 1 second.
3. Change to the low conductivity buffer and when the new conductivity level is stable, press  for 1 second.

The monitor is now ready for use

The setting of the 100%, and 0%, levels will remain in the transmitter until a new scale is set.

## 5.2 Menu

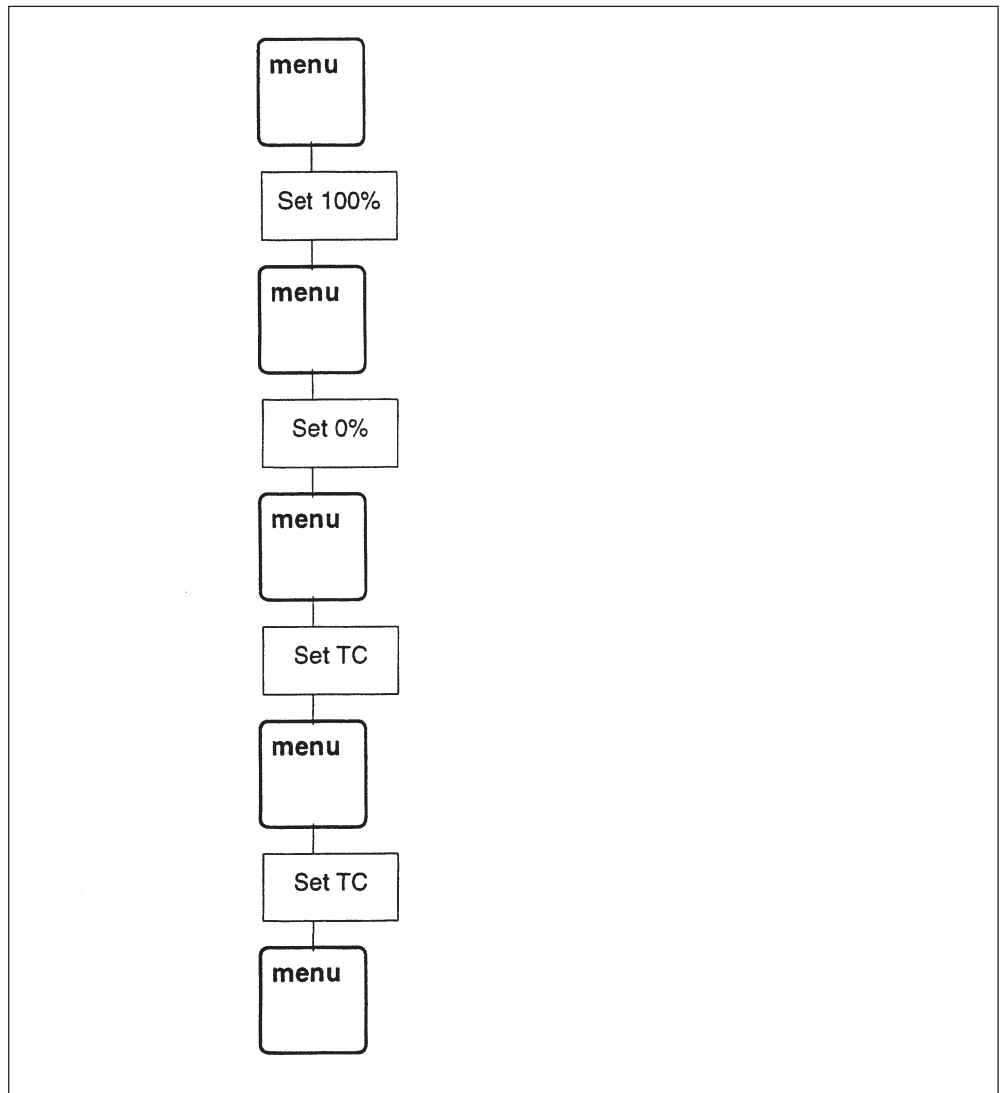

*Fig. 5 Menu overview.*

With the menu key, more advanced functions can be programmed, such as setting the conductivity to any range independent of the buffers and temperature compensation of the conductivity to reference temperature. The menu key is also used to step between the functions.

## 5.3 Fixed scale setting

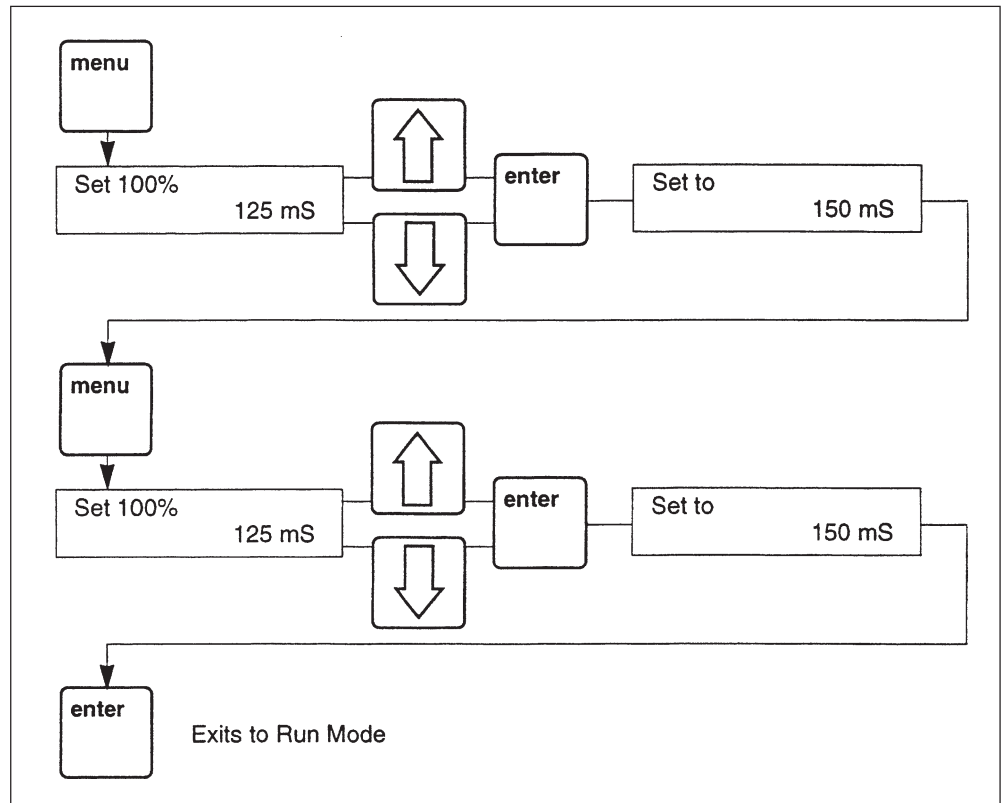











Fig. 6. Fixed scale setting

The 100% and 0% conductivity values can be set to any value of your choice, e.g. to extend the recorder full scale deflection and baseline to be slightly over and under the conductivity of the buffers. In addition it might be easier to read the conductivity values for specific peaks in a chromatogram if the scale is in tens or hundred of mS/cm.

In this example the conductivity scale is changed from a previously set 100% value of 125 mS/cm to a new value of 150 mS/cm. the 0% value is changed from 3.5 mS/cm to 0.00 mS.

1. Press , the previously set 100% value is displayed.
2. Use  and  keys to set the conductivity value to 150 mS.
3. Press  to accept the value. Set to 150 mS is displayed.
4. Press  again to set the 0% value.
5. Set the value for 0% to 0.00 mS using the  and .
6. Press  to accept the setting. Set to 0 mS is displayed.
7. Press  again to exit Run Mode.

## 5.4 Temperature compensation

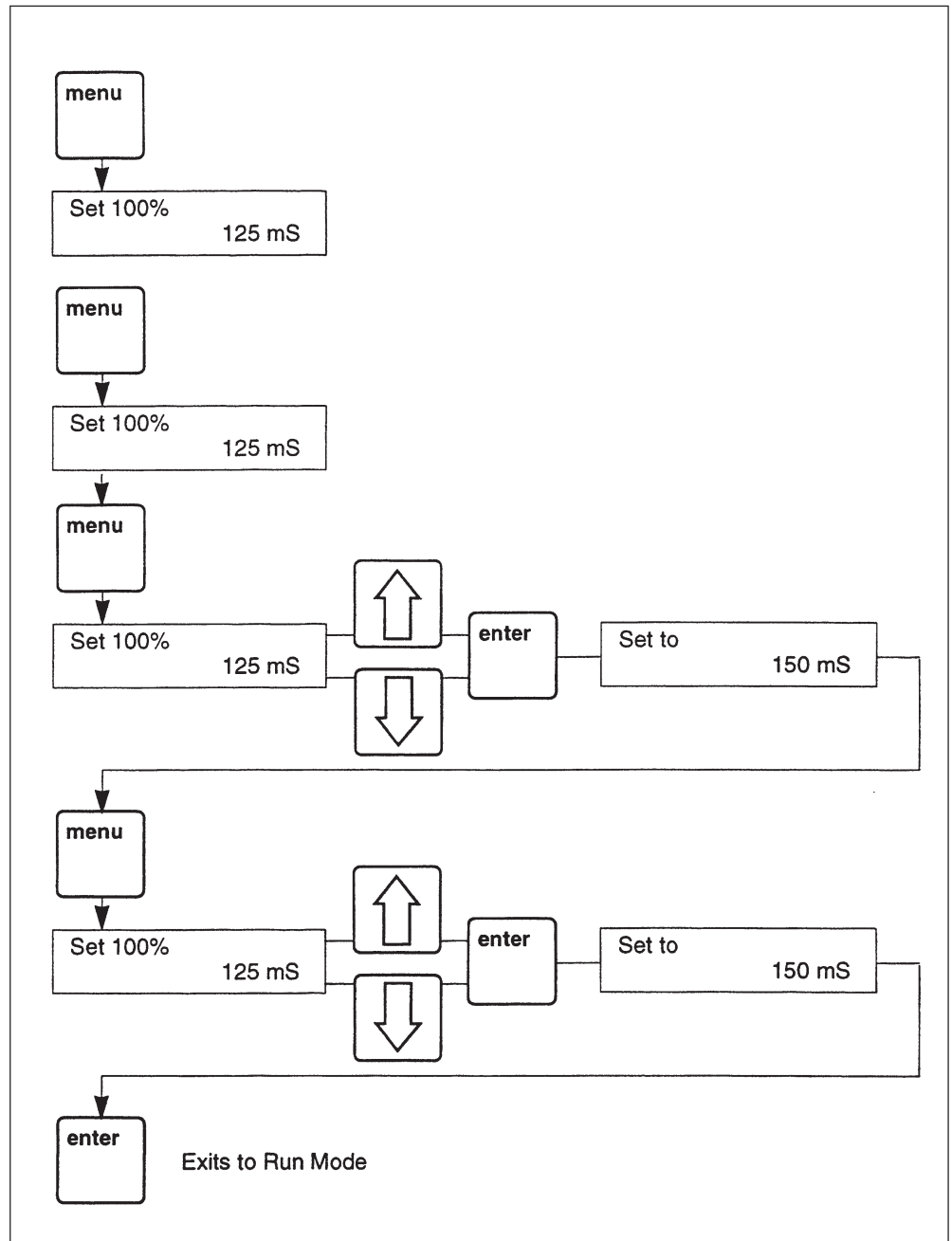











Fig. 7 Temperature compensation setting.

Temperature variations may influence conductivity. If high accuracy in the specific conductivity is required, a temperature compensation factor can be programmed together with a reference temperature. All conductivity values will then automatically be converted to the set reference temperature. If the temperature coefficient is unknown a recommended average temperature compensation factor of 2% (conductivity increase in percentage per °C) can be used for most common salts used in liquid chromatography applications.

When the temperature compensation is set to 0, this function is inactive.

In the example below the temperature compensation factor is changed from a previously set value of 2.1% to a new value of 2.0%. The reference temperature is change from 23.5 °C to 25.0 °C.

1. Press  three times until the Set TC menu is displayed.
2. Set the value of the temperature compensation factor to 2.0% using the  and  keys.
3. Press  to accept the value. Set to 2.0% is displayed.
4. Press  to set the temperature to which the conductivity values should be corrected. Tref 23.5 °C is displayed. Use the  and  keys to set the temperature to 25.0°C.
5. Press  to accept the temperature.
6. Press  again to exit to Run Mode.



## 5.5 Summary of commands in menu

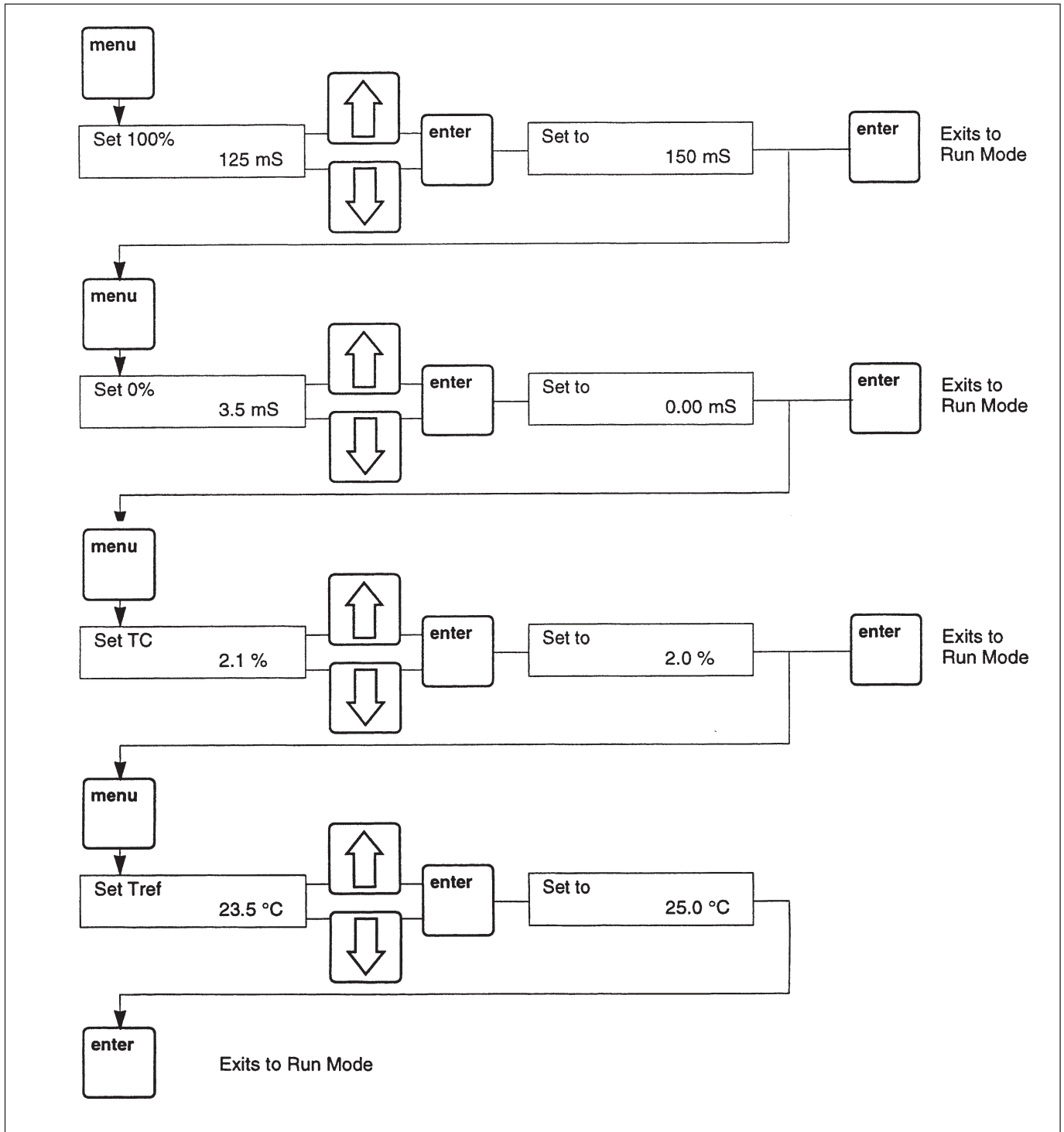


Fig. 8. Summary of commands.

## 6. Maintenance

The Conductivity Monitor requires very little maintenance if precautions are taken to protect the flow cell.

*Note:* When the flow cell is not in use, rinse the cell with de-ionised water, do not leave buffers to dry in the cell. When the monitor is stored for a longer period of time, fill it with 20% ethanol.

### 6.1 Cleaning the flow cell

The flow cell is designed to withstand cleaning-in-place with 1 M NaOH.

### 6.2 Calibration Mode

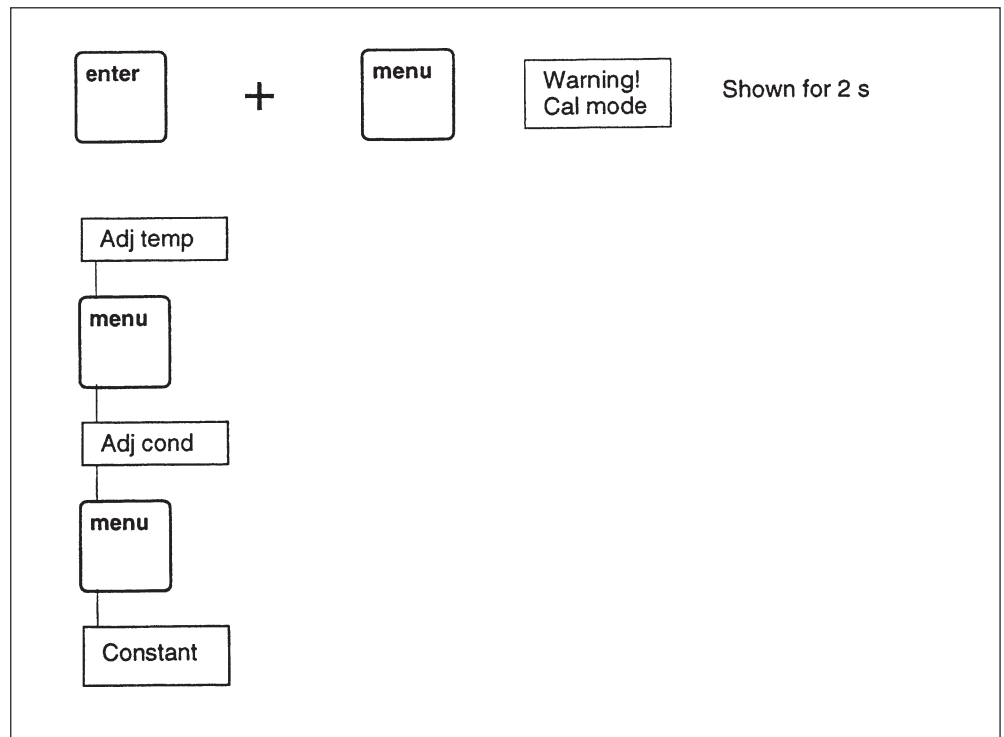


Fig. 9. Calibration Mode menu

## 6.3 Calibration of the temperature sensor

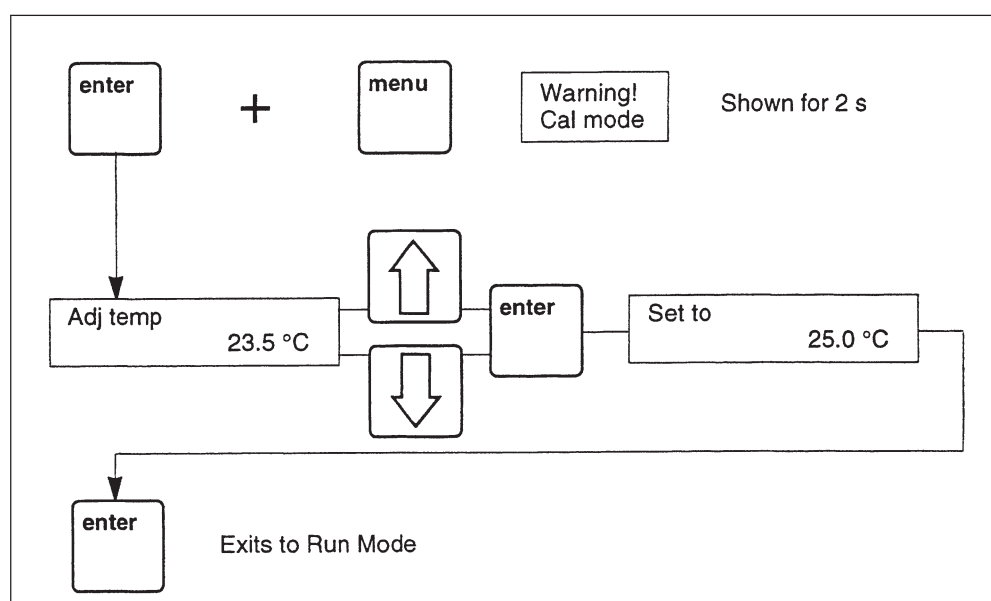


Fig. 10. Calibration of the temperature sensor.

The cell constant and the temperature sensor can be calibrated in the monitor calibration mode. New monitors are factory calibrated and it is not recommended to change these values.

Calibration of the temperature sensor is only necessary if the monitor is used to determine absolute conductivity with high accuracy.

In the example below the temperature is adjusted from 23.5 °C to the measured temperature of 25.0 °C.

1. Place the flow cell together with a precision thermometer inside a box or beaker to ensure that they are not exposed to draft. Leave them for 15 minutes to let the temperature stabilise.

2. Read the temperature on the thermometer.

3. Press  and  simultaneously for 2 seconds to enter the Calibration Mode.

4. The display shows Adjust temp. 23.5 °C. Use the  and  keys to set the temperature on the control unit to 25.0 °C the temperature shown on the thermometer.

5. Press  to accept the setting.

6. Press  again to exit the Calibration Mode.

## 6.4 Calibration of the cell constant

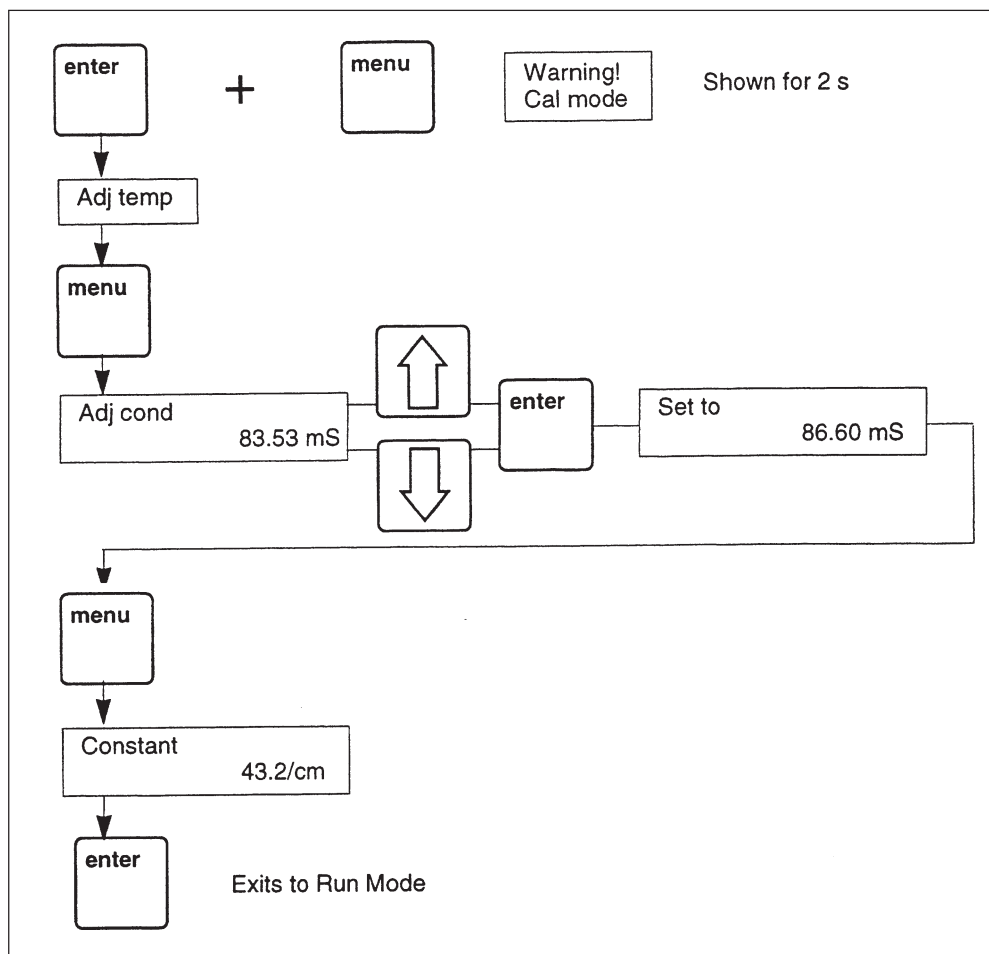


Fig. 11. Calibration of cell constant.

Normally it is not necessary to calibrate the flow cell as it is pre-calibrated on delivery. Calibration of the flow cell is only necessary when the monitor is to be used to determine specific conductivity with high accuracy.

In the example below the cell constant is changed by setting the conductivity value of a calibration buffer to theoretical value at an exact temperature.

1. Prepare a calibration solution of 1.00 M NaCl, 58.44 g/l.
2. Fill the flow cell with the calibration solution, either via a syringe or by pumping at least 15 ml through the cell.
3. Stop the flow and wait for 15 minutes, until the temperature is constant in the range 20-30 °C.
4. Read the conductivity value displayed.

*Note:* Make sure the Temperature Compensation is off.

5. Press  and  simultaneously for 2 seconds to enter the Calibration Mode.

6. Read the temperature in the Adj. temp menu. Compare the earlier displayed conductivity value with the theoretical value from the graph below, at the temperature of the calibration solution. If the theoretical and the displayed value correspond no further action is required. If the values differ proceed as outlined below.

7. Press  again until the display shows adjust cond. Use the  and  keys to enter the theoretical value.
8. Press  to accept the value. The new cell constant is automatically calculated.
9. Press  again and the new cell constant value will be displayed in the me  
Constant.  
*Note:* Do not change this value.
10. Press  to exit to run mode.

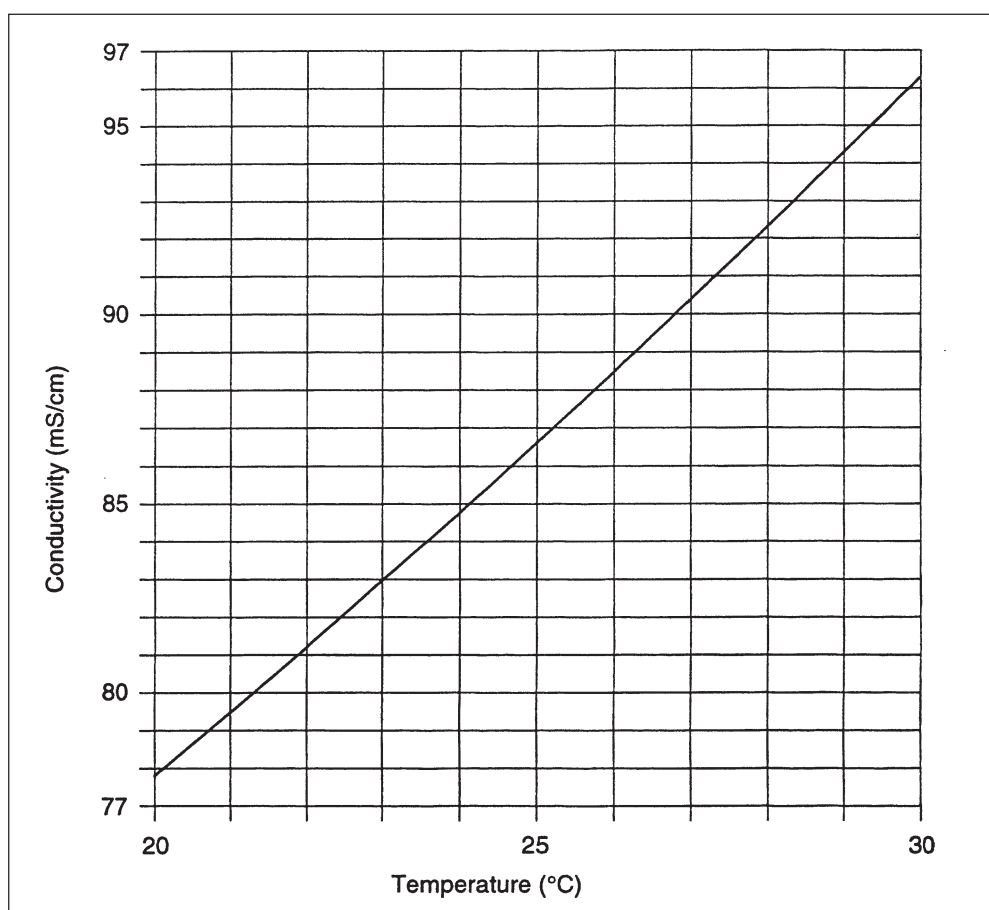


Fig. 12. Conductivity of 1 M NaCl at 20-30 °C.

## 6.5 Replacing the flow cell

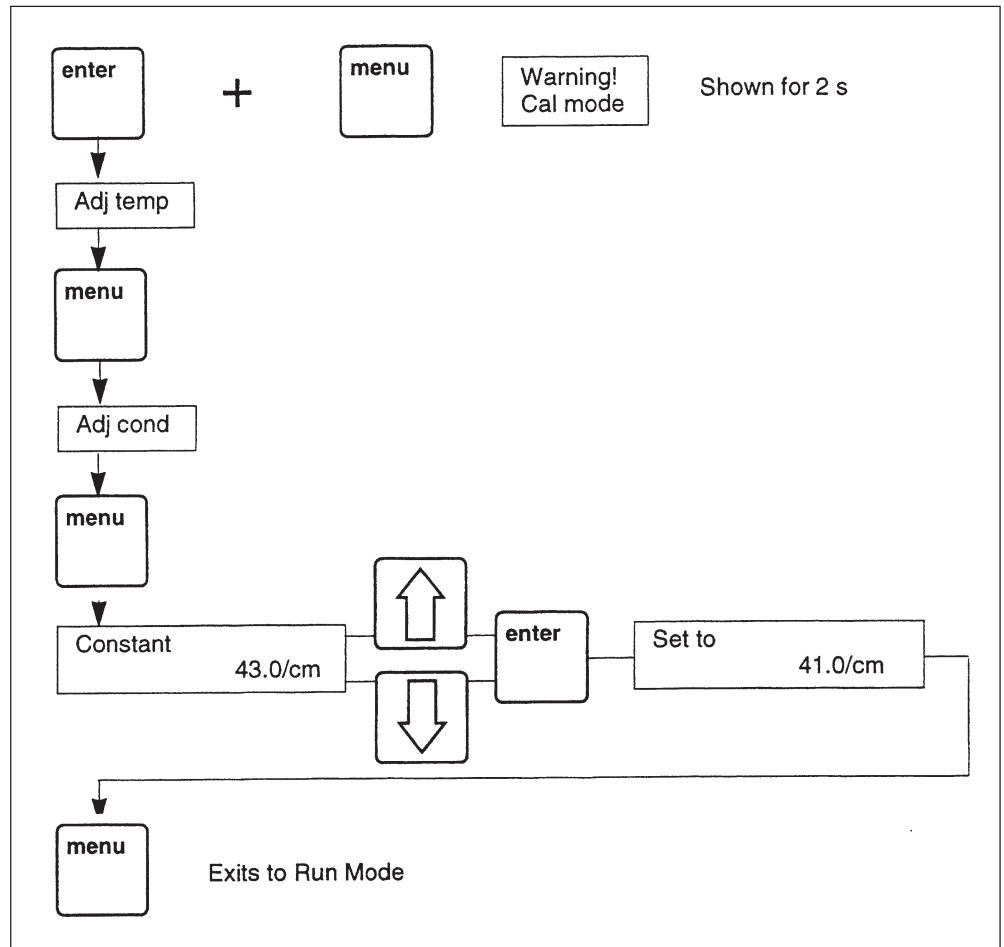









Fig. 13. Changing the cell constant.

If the flow cell is replaced with a new flow cell, the monitor must be calibrated with the new cell constant value. The new cell constant is written on the package. In the example below the cell constant is changed from 43.0/cm to 41.0/cm.

1. Press  and  simultaneously for 2 seconds to enter the Calibration Mode.
2. Press  twice until Constant 43.0/cm is shown in the display. Use the  and  keys to set the cell constant to 41.0/cm.
3. Press  to accept the value.
4. Press  again to exit to Run Mode.

## 6.6 Summary of commands in Calibration Mode

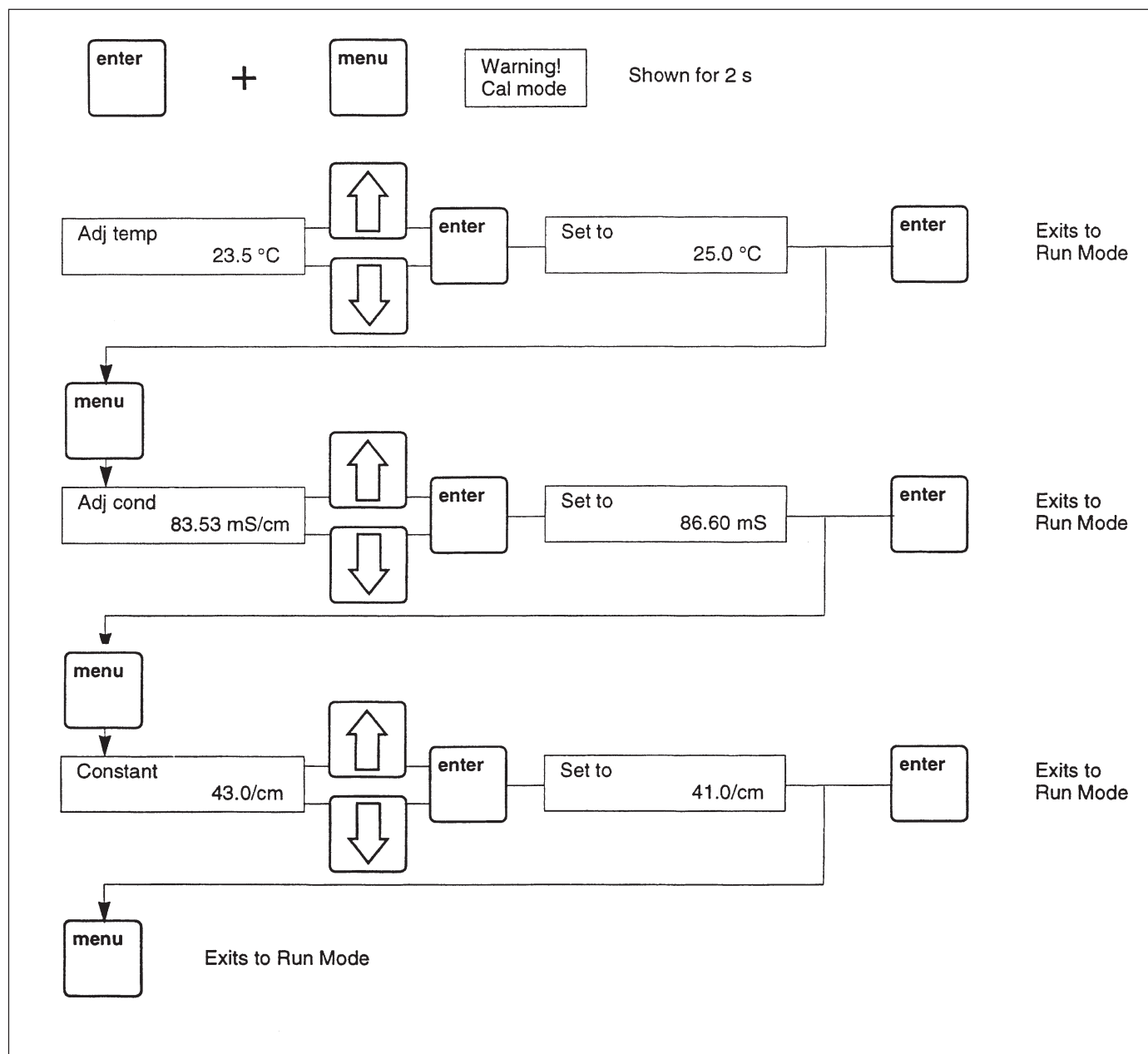


Fig. 14. Summary of commands.

# 7. Trouble Shooting

<b>Problem</b>	<b>Possible cause</b>	<b>Remedy</b>
Baseline drift or noisy signal monitor	Air in flow cell	Use a flow restrictor after the conductivity flow cell
described in	Dirty flow cell	Check for leaking tubing connections Clean flow cell according to procedure Maintenance 6.1
	Column not equilibrated or dirty	Equilibrate and/or regenerate column
The absolute conductivity value is wrong	The electrical grounding of the electrode is influenced by another instrument	Turn the flow cell so the end with the screws is electrode is influenced by facing the other instrument
	Temperature compensation is on	Set the temperature compensation factor to 0
Ghost peaks appear in the gradient profile	Charged sample detected (e.g. protein)	None required
	Air bubbles passing through flow cell	Use a flow restrictor after the conductivity monitor flow cell Check for loose tubing connections
Conductivity measurements with the same buffer appear to	Dirty flow cell	Clean flow cell according to procedure described in section 6.1
	Room temperature has decreased	Use temperature compensation factor
The monitor does not respond to conductivity change	Faulty flow cell or transmitter	Call a Amersham Biosciences service representative



## 8. Accessories

Designation	Code No.
Conductivity Flow Cell	19-1515-15

### Trademarks

Amersham and Amersham Biosciences are trademarks of Amersham plc.

© Amersham Biosciences AB 2000 – All rights reserved

All goods and services are sold subject to the terms and conditions of sale of the company within the Amersham Biosciences group which supplies them.

A copy of these terms and conditions is available on request.

#### **Amersham Biosciences AB**

Björkgatan 30, SE-751 84 Uppsala  
Sweden

#### **Amersham Biosciences UK Limited**

Amersham Place, Little Chalfont  
Buckinghamshire, HP7 9NA  
England

#### **Amersham Biosciences Inc**

800 Centennial Avenue, P.O. Box 1327  
Piscataway, New Jersey 08855  
USA

#### **Amersham Biosciences Europe GmbH**

Postfach 5480  
D-79021 Freiburg  
Germany



Amersham  
Biosciences



