

FOLIGAIN.L12x[®]

Professional Laser Comb

QUICK START GUIDE



Powered by: **HAIRMAX[®]**

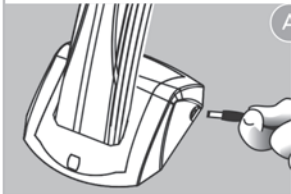
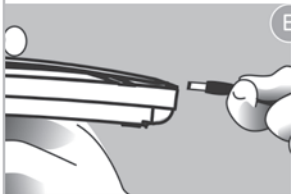
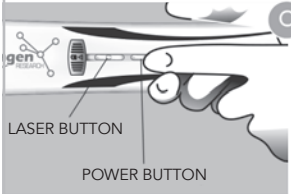
The FOLIGAIN.L12x® Professional Laser Comb powered by HairMax® is a very effective device when used correctly. It is manufactured using precision components and should last 10-15 years under normal usage. **Please take time to read the entire User Manual before operating the Laser Comb.**

FOLIGAIN.L12x® Professional Laser Comb



Use the Laser Comb for 8 minutes, 3 times a week. The treatments should be spaced every other day, for example, on Monday, Wednesday and Friday. Move the Laser Comb slowly 1/2 inch (1 cm) every 4 seconds.

It is recommended, but not essential, that the Laser Comb be used on a clean scalp after washing. Your hair can be wet, towel-dried or dry. The Laser Comb should be used prior to applying any topical treatments, gels, sprays or concealing products.

  <p>Figure A shows the Laser Comb being placed into a charging cradle. Figure B shows a hand inserting a power cord into the handle end of the Laser Comb.</p>	<p>Charge your Laser Comb using either the charging cradle (Figure A) or by placing the power cord directly into the handle end of the Laser Comb (Figure B). For initial charging, allow 2-3 hours to fully charge your Laser Comb. The battery indicator light of the Laser Comb will be yellow and will turn green when fully charged. When the Laser Comb requires charging, the battery indicator light will flash red and the battery icon will appear on the information display.</p>
 <p>Figure C shows a hand pressing the power button and the laser button on the Laser Comb. Labels indicate the "LASER BUTTON" and "POWER BUTTON".</p>	<p>To turn on the Laser Comb, press the power button and hold until power button light turns green. Next, press the laser button, which will light up red. That means the Laser Comb is ready for use (Figure C).</p>

Laser Light - Avoid Direct Eye Exposure

The FOLIGAIN.L12x® Professional Laser Comb is indicated to treat Androgenetic Alopecia, and promote hair growth in males who have Norwood Hamilton Classifications of IIa to V and in females who have Ludwig (Savin) I-4, II-1, II-2 or frontal patterns of hair loss and who both have Fitzpatrick Skin Types I to IV.



Before beginning treatment, be sure to test the Laser Comb for sensitivity (see page 5 of your User Manual). To start your treatment, place the Laser Comb flat on your scalp so that both rows of teeth are touching the scalp (Figure D).



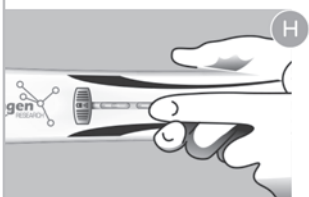
Start at your hairline. Move the Laser Comb from front to back of scalp in a slow combing motion, stopping every $\frac{1}{2}$ inch (1 cm) and pausing for 4 seconds (Figure E). Do not lift the Laser Comb off your scalp while combing from front to back.



Next, lift the Laser Comb off your scalp and start at a new position at the hairline. Continue in this manner from front to back, then side to side, (Figure F) until you cover your entire scalp with multiple passes. When the Laser Comb is on, it will beep and vibrate every 4 seconds. This means that it is time to move the Laser Comb to the next area of your scalp. Repeat these passes for the total treatment of 8 min.



If you have long hair, use your fingers to lift your hair out of the way while using the Laser Comb to be sure that the laser energy reaches your scalp (Figure G).



After treatment, turn off the Laser Comb by pressing and holding the power button (lower button) until the green indicator light turns off (Figure H). Store the device in the charging cradle for future use. You may leave the power cord plugged into the charging cradle at all times.

PLEASE NOTE:

- The total time required for treatment with the FOLIGAIN.L12x® Professional Laser Comb is 8 minutes.
- It is important to move the Laser Comb very slowly. Remember, you are bathing the scalp with laser light. If you move it too fast, the treatment will not achieve the optimal results to promote hair growth.
- While there are no negative effects from overuse, it is recommended that you use the Laser Comb according to directions.
- The Laser Comb has advanced electronics to make sure that the device is working properly. If you are experiencing difficulties with your device, please refer to the Troubleshooting section of your User Manual before sending it back to us for repair.

FOR INTERNATIONAL TRANSLATIONS VISIT:

www.anagenresearch.com

- FR Pour une traduction internationale, visitez: www.anagenresearch.com
- ES Para traducciones internacionales, por favor visite: www.anagenresearch.com
- IT Per traduzioni internazionali, per favor visitate: www.anagenresearch.com
- DE Für internationale Übersetzungen besuchen Sie bitte: www.anagenresearch.com
- PT Para as traduções internacionais, visite: www.anagenresearch.com

REGISTER

To validate your warranty
please register online at
www.hairmax.com/warranty

CONTACT

Anagen Research, LLC
8560 West Sunset Boulevard, Suite 500
West Hollywood, California 90069
www.anagenresearch.com

