



Philips Domestic Appliances and Personal Care

Service Manual

PRODUCT INFORMATION

Small-volume brush head	: V-shape following the contours of the teeth.
Active tip	: Independently moving brush with longer and more flexible bristles.
Controlled Pressure System	: Brush flexes back when too much pressure is applied.
End-rounded bristles	: Gentle to gums and teeth.
2 speeds	: The second speed is for more gentle cleaning in sensitive areas.
Duo-timer	: The first 2 minutes of the brushing process, the green pilot light of the Duo-timer blinks. After these two minutes, the brushing process is briefly interrupted and the green pilot light burns continuously. This interruption will occur every two minutes. The Duo-timer measures the actual brushing time, even if the appliance is switched off during the brushing process for max. 10 seconds.
Almost Empty Indication	: If the battery is almost empty, the pilot light is red and blinks for 10 seconds after the toothbrush has been switched off.
Full Charge Indication	: If the battery in the toothbrush is being charged, the pilot light is green and blinks. If the battery has been fully charged, the pilot light burns continuously
Quick charge	

TECHNICAL DATA

Voltage	: 200 - 240 V, 50 - 60 Hz
Power consumption	: 3 Watts
Speed (Top)	: up to 15.000 brush strokes per minute
Pressure control click at	: 250 gr.
Charging time	: empty to full: 3 hours
Operating time	: full to empty: 25 minutes
Battery type	: NiCd

SERVICE INFORMATION

Hand unit tooth brush (1)

Never open the hand unit of the tooth brush, because the unit is not water tight anymore.
In case of a defect, the complete unit must be replaced.

Charging unit (2)


In case of a defective charging unit or cordset, the complete charging unit must be replaced.

NOTES

When discarding the appliance, remove the built-in rechargeable battery.
The battery should be completely empty of charge.

Product meets the requirements regarding interference suppression on radio and TV.

EXPLODED VIEW + PARTS LIST HX 2555

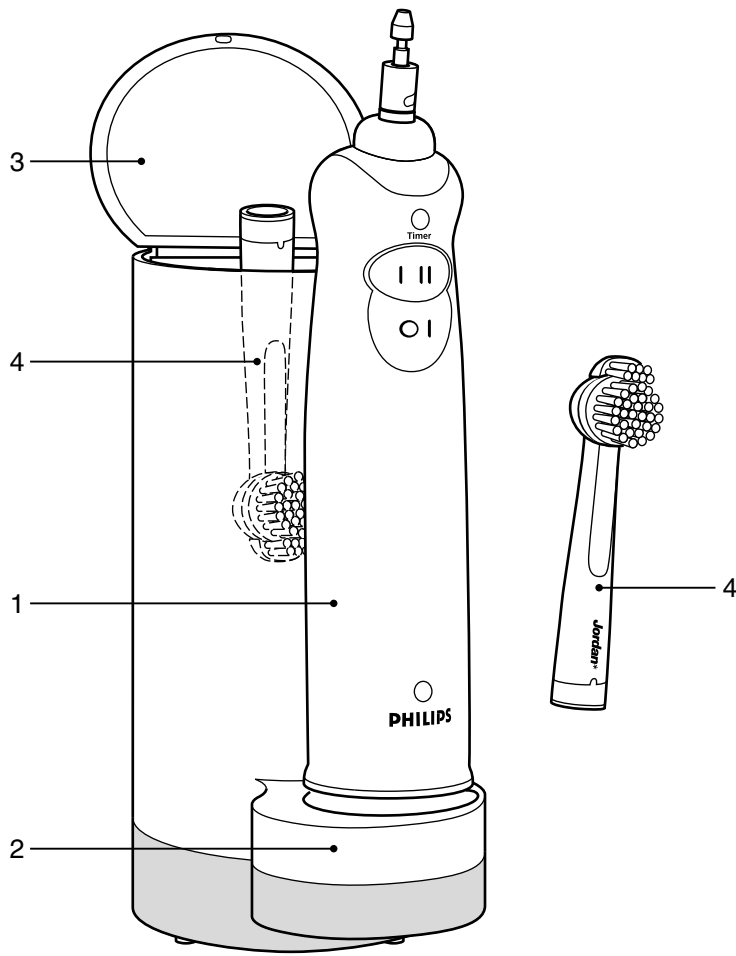

Pos	Description	Service code
1	Rechargeable hand unit (NiCd)	4203 035 71600
2	Charging unit 100-240 V, EURO plug	4203 035 71620
3	Cover Container	4822 442 01902
4	Set of brushes, 2 pcs in blister pack	4822 479 20272
	Dir. for use: UK, D, F, NL, E, P and TR	4203 032 93530
	Dir. for use: UK, D, N, S, Dk, SF and Gr	4203 032 93580