

User's Manual

Digital Spirometer

IQmark

Document ID: PUB-116

Revision B



The names of any providers and patients used in illustrations or examples in this document are fictitious.

Every effort has been made to ensure this manual is accurate, complete, and useful. Please let us know if you have any suggestions for improvement using one of the following means of contact:

AFHCAN
Training Department
4000 Ambassador Drive
Anchorage, AK 99508

Phone: 877 885-5672
Fax: 907 729-2269
Email: customersupport@afhcan.org

Copyright © 2006 Alaska Native Tribal Health Consortium (ANTHC). All rights reserved. No part of this publication may be reproduced, stored in a retrieval system, or transmitted, in any form, or by any means, including photocopying, electronic, mechanical, recording or otherwise, without the prior written permission of Alaska Native Tribal Health Consortium.

Table of Contents

Section	Title	Page
Preface	v
	About this Document	v
	Related Documents	v
	For More Information	v
Section 1 – Introduction	1
1.1	Introduction to the Spirometer	1
1.1.1	Functional Description.....	1
1.1.2	Main Components.....	2
1.1.4	AFHCAN Software Features.....	4
1.1.4.1	How the AFHCAN Software Works with Midmark Software.....	4
1.1.4.2	Spirometry Screen	5
1.1.5	Features of the Midmark Spirometry Software	7
1.1.5.1	Introduction to Available Spirometry Tests	7
1.1.5.2	Obtaining Spirometry Data.....	11
1.1.5.3	Viewing Spirometry Data.....	16
1.1.5.4	Calibration Screen	19
1.1.5.5	Settings Dialog Box.....	23
1.2	Particulars of the AFHCAN Installation.....	27
1.2.1	Version Compatibility	27
1.3	Warnings and Cautions	27
Section 2 – Operation	29
2.1	Basic Operating Procedures	29
2.1.1	Starting the Spirometer	29
2.1.2	Calibration	29
2.1.3	Obtaining a Spirometry Report.....	31
2.1.3.1	Patient Preparation and Coaching	31
2.1.3.2	Conducting the Test.....	31
2.1.3.3	Reviewing and Selecting Tests.....	32
2.1.3.5	Saving the Spirometry Report to a Case.....	33
2.1.4	Viewing Spirometry Reports	33
Section 3 – Clinical Considerations	35
3.1	Guidelines for Clinical Success	35
3.2	Common Mistakes	35
3.3	Tips and Tricks	35
Section 4 – Routine Maintenance	37
4.1	Care and Cleaning.....	37
4.2	Replacement Items.....	37
4.3	Elementary Troubleshooting.....	37

List of Illustrations

Figure	Title	Page
1	Digital spirometer.....	2
2	Disposable mouthpiece and nose clip.....	3
3	Calibration syringe and adapter.....	3
4	Spirometry screen.....	5
5	Data entry screen for supplemental patient information	6
6	Introduction to three principal screens in Midmark software.....	7
7	The basic FVC curve.....	8
8	Data showing flow volume loop.....	9
9	Comparing a series of FVC tests.....	9
10	Example of VC test results.....	10
11	Example of MVV test results	10
12	Midmark data acquisition screen.....	11
13	Patient incentive displays.....	15
14	Examples of patient feedback messages.....	16
15	Midmark data display screen – numerical data	17
16	Midmark data display screen – graphical data	17
17	Midmark calibration screen – first calibration	20
18	Midmark calibration screen – subsequent calibrations.....	21
19	Calibration pump window – single flow	22
20	Calibration pump window – multiflow test underway	22
21	Available tabs in settings dialog box.....	23
22	Settings tabs – first row	24
23	Settings tabs – second row	25
24	Settings tabs – third row.....	26
25	Double-teaming the calibration syringe	30

List of Tables

Table	Title	Page
1	Spirometer screen buttons and their functions	6
2	Items in the test selection and control panel.....	11
3	Midmark screen function buttons and applicable screens.....	13
4	Data display screen items that are active when a "spirometry session" is open.....	18
5	Troubleshooting guidelines	37

Preface

About this Document

This document is part of the set of user manuals provided with the AFHCAN Cart. These user manuals, covering various topics, are normally assembled into a binder delivered with each Cart. This modular design has the following advantages:

- the set of manuals provided with your cart includes documents for the specific peripheral devices installed
- each document is a stand-alone publication, so as new devices or features are added to the Cart, new manuals can be added to the existing binder
- user information that is common to all items of equipment does not need to be repeated in each module, but can be covered in separate modules and referenced as needed

Related Documents

This document assumes you have read the introductory hardware and software manuals included in this binder.

The original manuals provided with the equipment were included in a set of materials delivered with the AFHCAN Cart. Those manuals can be used to supplement the information provided in this document. Be aware, however, that items installed on an AFHCAN Cart may have been modified slightly, so the features as described in the original product manuals may not apply.

For More Information

This document describes the equipment to a level of detail that will meet most user's needs in the context of clinical use of the AFHCAN Cart. For more information, contact AFHCAN Customer Support:

AFHCAN Customer Support
Phone: 888 449-4435
Fax: 907 729-2269
email: customersupport@afhcan.org

User manuals, specifications, brochures, and other information that may be helpful can be obtained at Midmark's website:

www.midmark.com

Section 1 – Introduction

1.1 Introduction to the Digital Spirometer

1.1.1 Functional Description

Spirometry measures the patient's ability to move air into and out of the lungs.

Spirometry is a medical screen test that measures various aspects of breathing and lung function. It is performed by using a spirometer, a special device that registers the amount of air a subject inhales or exhales and the rate at which the air is moved into or out of the lungs. Spirograms are tracings or recordings of the information obtained from the test. The most common spirometric tests require that the subject exhale as forcefully as possible after taking in a full, deep breath. The subject's effort is called the forced expiratory maneuver.¹

The IQmark Digital Spirometer, working with Midmark software, measures the flow rate of air exhaled through a disposable mouthpiece attached to a hand-held sensor device. The Midmark software tracks the flow rate over time and calculates total volume of air expelled. The data obtained is displayed both graphically and numerically.

The most common spirometric test is the forced vital capacity (FVC) test. As described above, the patient inhales as deeply and fully as possible, and exhales through the mouthpiece as vigorously and completely as possible. Although the FVC test is most common, additional types of tests can be performed (refer to Midmark's Operation/Service Manual supplied with the device for details). In addition to recording flow rates and volumes for exhalations, the digital spirometer can record the same data for inhalations.

Normal flow rates and volumes have been established scientifically. Most healthy patients can expel 70 to 80 percent of their lung capacity in the first second. The Midmark software obtains predicted performance criteria based on the patient's age, height, weight, gender, and race, and then compares the patient's actual performance against these criteria. In addition, the Midmark software gauges the quality of the patient's effort, and provides messages on how to improve performance.

The Midmark software and the AFHCAN software work together. The Midmark software controls the digital spirometer, and captures and displays the data. The AFHCAN software activates the Midmark software, provides some patient information, and associates the spirometry data with a case. Once the Midmark data is in a case, it becomes a permanent record which can be sent to another provider for review, printed out for the patient's records, and eventually archived electronically.

It has been shown that good results depend on the patients making their best possible effort. Patients must be cooperative, comfortable, and very clear on how to perform the maneuver.

¹ NIOSH Spirometry Training Guide, December 12, 2003.
www.cdc.gov/niosh/docs/2004-154c/pdfs/2004-154c-ch2.pdf

The Midmark software can generate a tentative diagnosis based on how the spirogram compares with the predicted norms for the patient. The automatic diagnosis is a general guide that must be confirmed or modified by the physician responsible for the patient. Please read and understand the following statement, which is from the Midmark Operation/Service Manual:

Physician's Responsibility

The IQmark Digital Spirometer can provide test interpretations if the user enables the Auto Interpretation feature. These interpretations are for the exclusive use of licensed physicians or personnel under their direct supervision. The suggested interpretation and the numerical and graphical results should be examined with respect to the patient's overall clinical condition. Final analysis should always be determined and verified by a physician.

Spirometry is an effort dependent test. Proper administration of the test is the physician's responsibility, as is making a diagnosis, obtaining expert opinions on the results, and implementing the correct treatment, if indicated.²

1.1.2 Main Components

The spirometer consists of the following components (Figures 1 through 3):

- spirometer handle with sensing electronics and computer cable
- mouthpiece (called the *disposable pneumotach mouthpiece* or DPM by Midmark)
- disposable nose clip (optional, but highly recommended for good results)
- calibration syringe (with adapter)

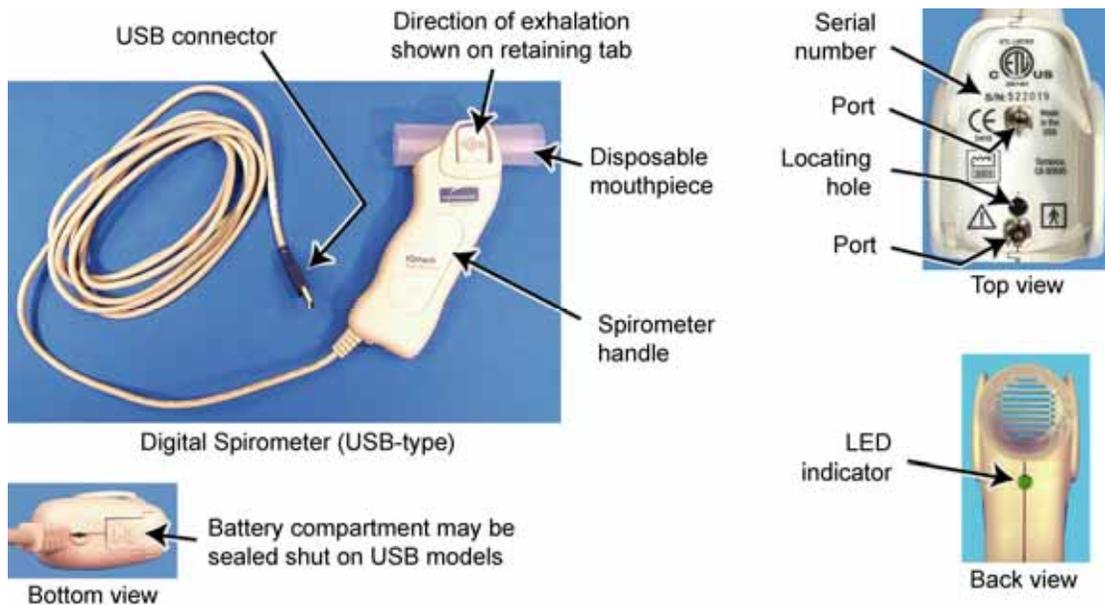


Figure 1
Digital spirometer

² Midmark Diagnostics Group. IQmark Digital Spirometer Operation/Service Manual, Version 7.0. p iii

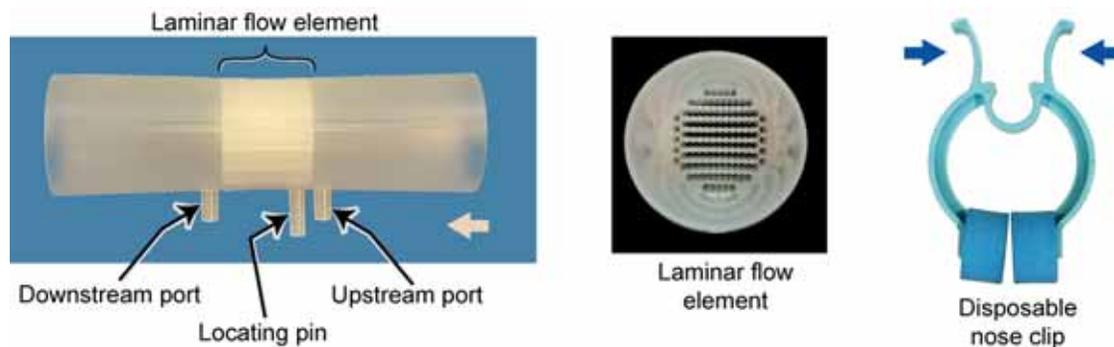


Figure 2
Disposable mouthpiece and nose clip

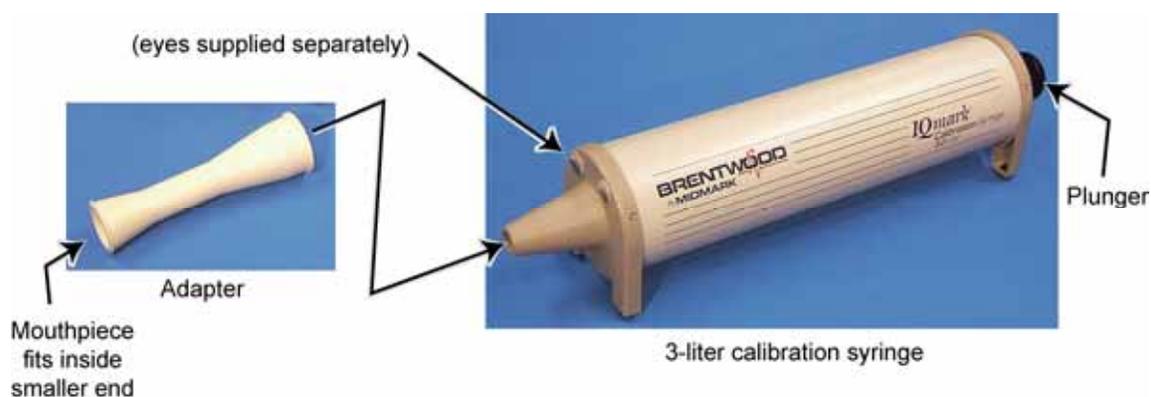


Figure 3
Calibration syringe and adapter

The spirometer handle is small, lightweight, and easily manipulated by the patient. At the top is a slot which receives the mouthpiece and secures it in place. On the back of the handle is an LED which illuminates green when the device is communicating with the computer.

The mouthpiece is a low-cost, single-use, disposable item. Inside the mouthpiece there is a laminar flow element (the corrugations in the middle of the mouthpiece) and two pressure ports, one on either side of the laminar flow element. The laminar flow element performs three functions:

- it offers a small amount of resistance to air flow, and this resistance to air flow creates a small difference in air pressure between the upstream port and the downstream port
- it smoothes or stabilizes the air flow through the mouthpiece so good pressure readings can be obtained
- it allows air to be exhaled with virtually no perceived resistance as far as the patient is concerned

The spirometer handle actually measures the difference in air pressure between the upstream port and the downstream port. The amount of pressure difference is related to the amount of air being forced through the mouthpiece at any given time. The air upstream of the laminar flow element will be at a slightly higher pressure than the air downstream from it. By monitoring the pressure

difference and tracking it over time, the Midmark software can compute flow rates and flow volumes.

The Midmark calibration syringe is used to verify instrument performance.

The mouthpiece includes a locating pin next to the upstream pressure port which ensures that the mouthpiece can only be installed one way. Ridges on the sides of the mouthpiece fit between the retaining tabs on the spirometer to ensure correct alignment and complete insertion into the pressure ports. Arrows on the tabs indicate the direction that exhaled air will travel.

1.1.4 AFHCAN Software Features

1.1.4.1 How the AFHCAN Software Works with Midmark Software

The AFHCAN software interacts with the Midmark software in three ways:

- activates the Midmark software
 - when creating a spirometry report and adding it to a case
 - when viewing a spirometry report already saved into a case
- feeds information from a case to the Midmark software
- configures the Midmark software so that some features are active and some are disabled

The AFHCAN software provides the following information to the Midmark software.

- patient's first name, middle initial, and last name
- patient's social security number
- patient's date of birth
- patient's gender
- provider's last name and first name (the user currently logged into AFHCAN software)

Midmark designed its software to be interactive with other software packages (such as the AFHCAN software) and therefore gave programmers the ability to activate or deactivate certain features of their software. A number of screens and features described in Midmark's Operation / Service Manual do not apply because the AFHCAN software handles those functions (for example, archiving a case). In addition, some features are disabled to ensure that the data in a case, which is an official medical record, is not accidentally altered or deleted at a later time. The AFHCAN software enables and disables features in the Midmark software as needed to ensure the integrity and security of medical records.

1.1.4.2 Spirometer Screen

Pressing the **Spirometer** button on the *Add To Case* screen brings up the *Spirometer* screen (Figure 4).

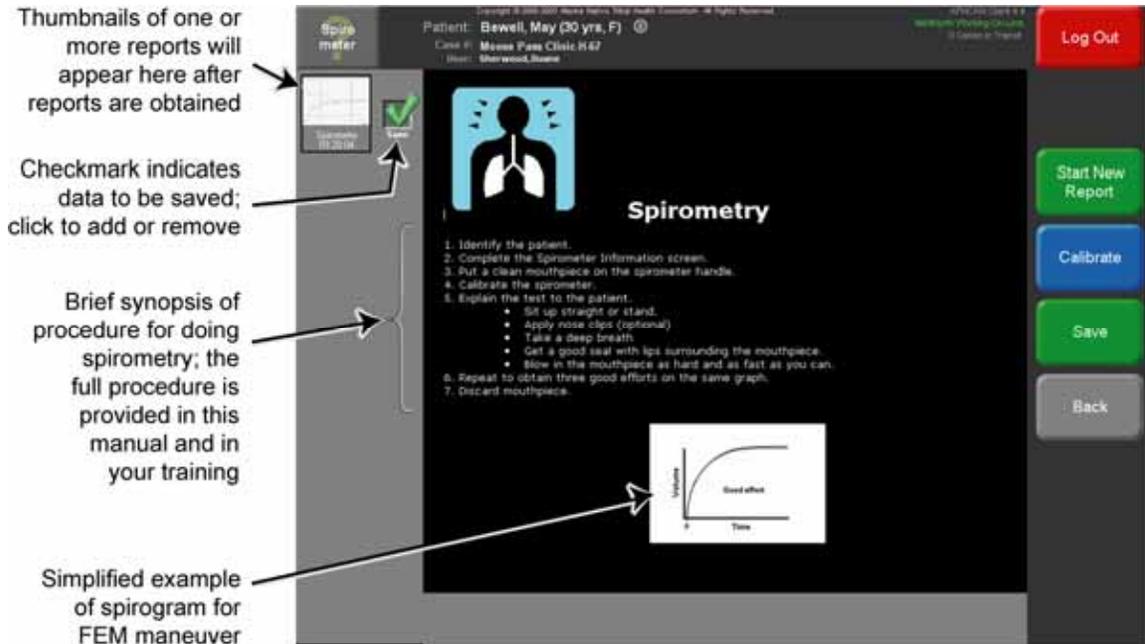


Figure 4
Spirometer screen

Accessing the *Spirometer* screen can be thought of as opening a spirometry "session." The session remains open as long as the case is open. The case remains open until it is sent, placed on hold, archived, or discarded. During the session, one or more spirometry reports can be created. Furthermore, an individual spirometry report can be accessed from the *Case* screen, and more spirometry tests can be added to that report.

A typical sequence will involve pressing the **Start New Report** button. The AFHCAN software automatically provides the Midmark software with patient identification information and some of the information needed for predicting performance. However, additional information must be supplied. Because the information such as weight, height, and smoking habits can change, the information must be entered each time a spirometry session is initiated (in other words, each time the **Spirometer** button is pushed on the *Add To Case* screen). When the **Start New Report** button is pressed, the screen shown in Figure 5 comes up.

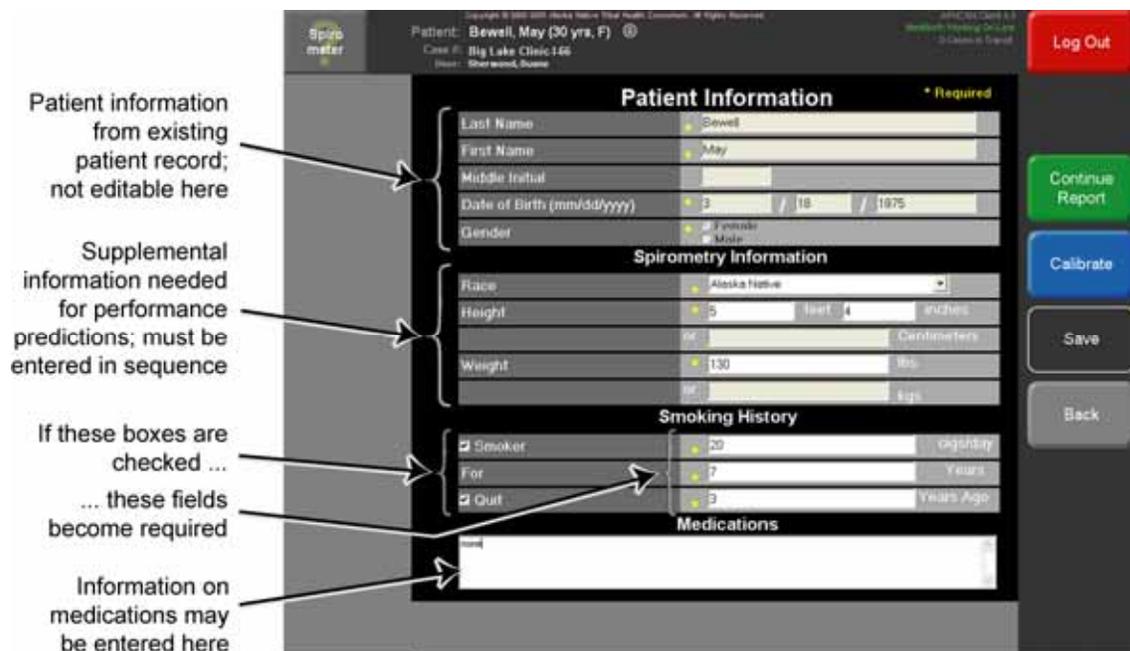


Figure 5
Data entry screen for supplemental patient information

Supplemental patient information has to be entered only once while a case is open. Once added, it is automatically incorporated into all the reports created in that session. On completion of the data entry, pressing **Continue Report** proceeds directly to the spirometry test screen in the Midmark software. Table 1 shows the buttons that appear on the *Spirometer* screen and their functions.

Table 1
Spirometer screen buttons and their functions

Button	Function
	If supplemental patient information has not yet been added in this session, the Start New Report button brings up the data entry screen for supplemental patient information. If supplemental patient information has already been added in this session, pressing the Start New Report button brings up the Midmark spirometry software.
	This button replaces the Start New Report button when the data entry screen for supplemental patient information is displayed. Pressing Continue Report stores the supplemental information and opens the Midmark spirometry software.
	Brings up the Midmark calibration screen. It is good practice to calibrate on a regular basis to verify equipment operation. Daily calibration is recommended by the American Thoracic Society.

Button	Function
	Saves all spirometry data with a green checkmark into the case, and returns to the Case screen. If no spirometry data has been obtained when the Save button is pushed, a message to that effect will be displayed and the software will return to the Spirometer screen.
	Returns to the Add To Case screen. If there are spirometry thumbnails in the left bar area when the Back button is pressed, a message will advise that data will be lost.

1.1.5 Features of the Midmark Spirometry Software

The Midmark software has three principal screens (Figure 6):

- data acquisition screen
- data display screen
- calibration screen

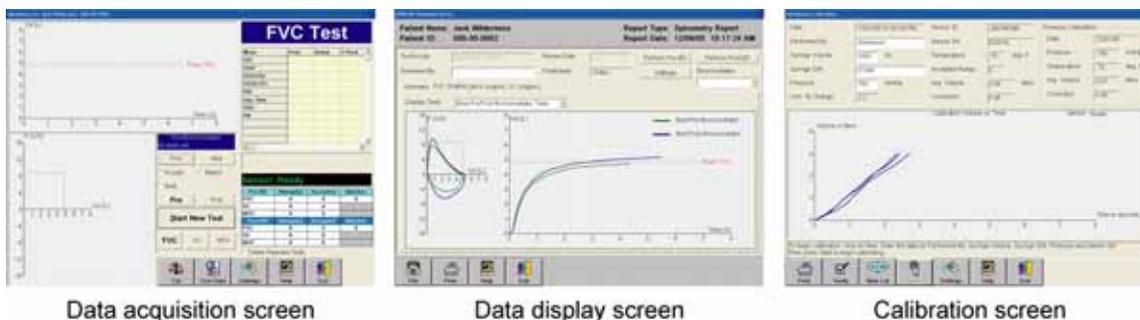


Figure 6

Introduction to three principal screens in Midmark software

The following functions are provided by or are accessible in the Midmark software:

- analyzing and capturing one or more spirometry reports
- printing a spirometry report
- reviewing a spirometry report (works with AFHCAN software to provide this feature)
- selecting various display and print formats
- changing some system settings (e.g., date format, English/metric)

1.1.5.1 Introduction to Available Spirometry Tests

The Midmark Digital Spirometer can perform a variety of tests and gather various items of data. Details on the use of different types of spirometry tests and data is beyond the scope of this manual. More information can be found in the Midmark Operation/Maintenance Manual, in the NIOSH Spirometry Training Guide (<http://www.cdc.gov/niosh/docs/2004-154c/>), or in other medical references.

This manual focuses on the most common spirometry test, the FVC or *forced vital capacity*. Although the FVC test is the primary focus of this manual, a basic familiarity with the nature of the other tests available will be helpful if those tests are requested.

Forced Vital Capacity Test

A single FVC test will produce two graphs (Figure 7):

- flow-volume curve (flow rate in liters per second vs. total volume expelled)
- volume-time curve (cumulative volume exhaled vs. time in seconds)

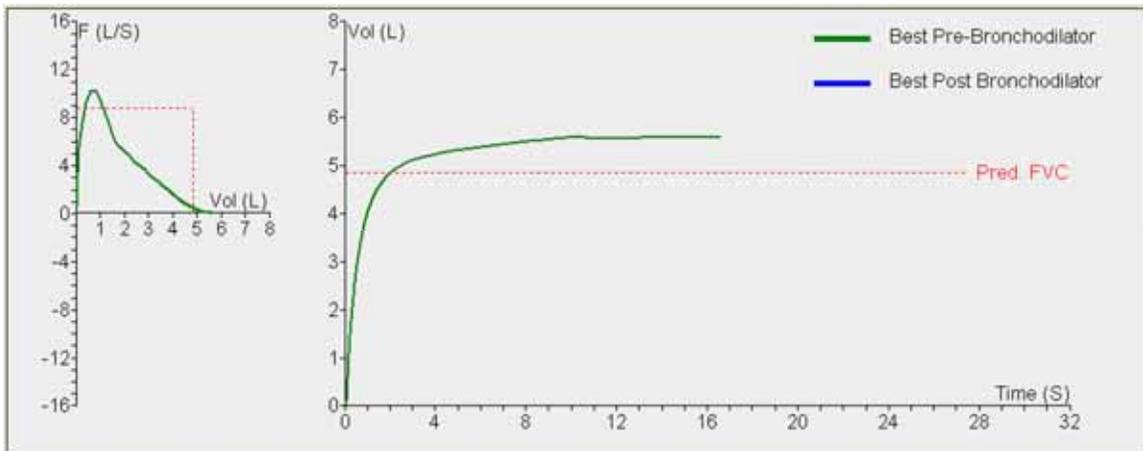


Figure 7

The basic FVC curve

The horizontal dotted line in the second curve is the predicted value of the forced vital capacity (FVC), which is the maximum amount of air that can be expelled by the lungs. In Figure 7, the patient exceeded the predicted value.

Recording Inhalation

If the patient inhales through the mouthpiece after the FVC test, inhalation data will be added to the graph, as shown in Figure 8. This information, called the *flow volume loop*, is shown in the negative part of the flow volume curve, and is automatically recorded.

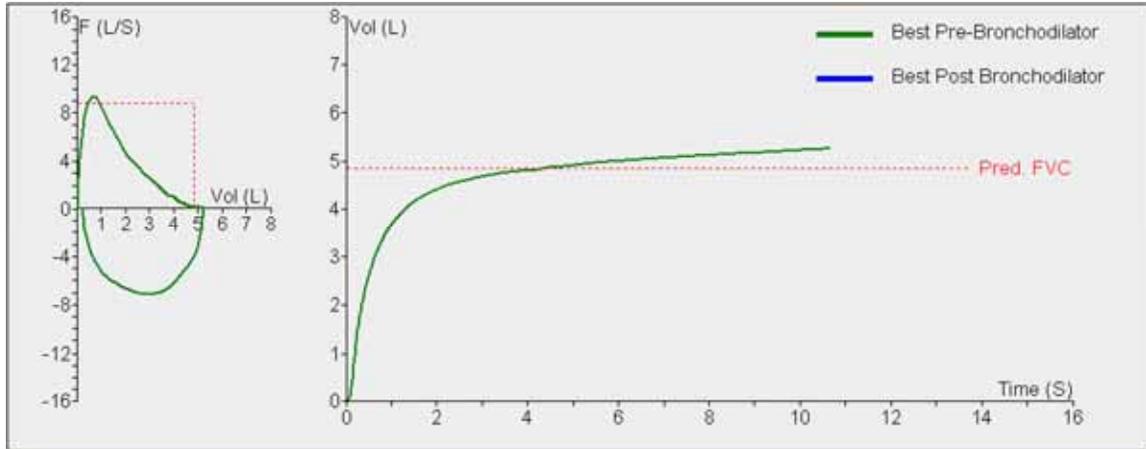


Figure 8
Data showing flow volume loop

Comparing Multiple Tests

A number of FVC tests (as well as other types of tests) can be added to a given spirometry report. Each spirometry report is stored as a separate thumbnail in an AFHCAN case. Multiple tests can be viewed on the same graph, as shown in Figure 9. If desired, the software allows tests to be identified as Pre Bronchodilator or Post Bronchodilator.

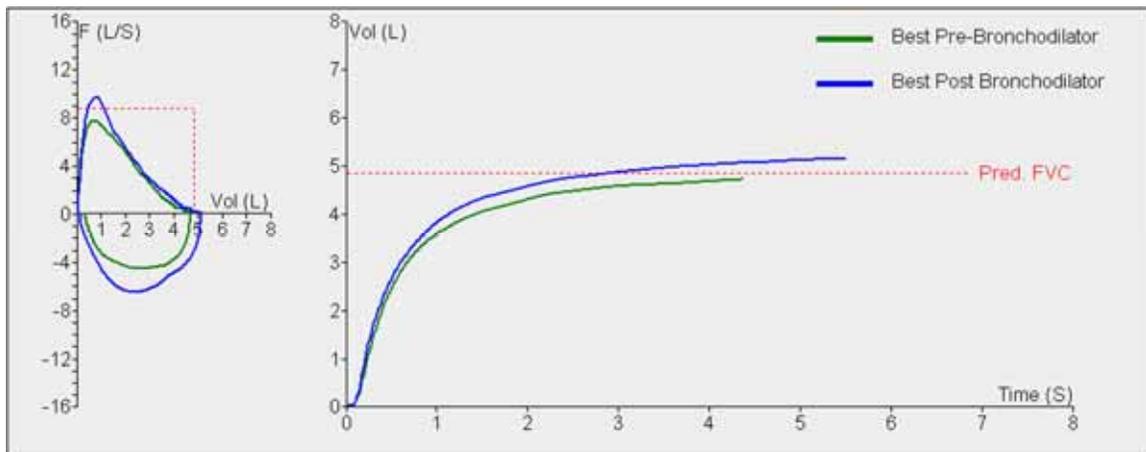


Figure 9
Comparing a series of FVC tests

Vital Capacity Test

The vital capacity (VC) test measures the full capacity of the lungs measured at slower rates of inhalation and exhalation. The test begins with the patient taking a full breath, exhaling into the mouthpiece, a few tidal breaths (normal breaths), and finally a full inhalation and full exhalation. The final inhalation and exhalation should not be done too slowly, but neither should it be forced. If the test is completed in about 25 seconds or less, all the breaths will be displayed on the graph, as shown in Figure 10.

The VC test is not normally required in the clinical situations where the AFHCAN Cart is usually installed. Additional details on performing this test are included in the Midmark Operation / Service Manual.

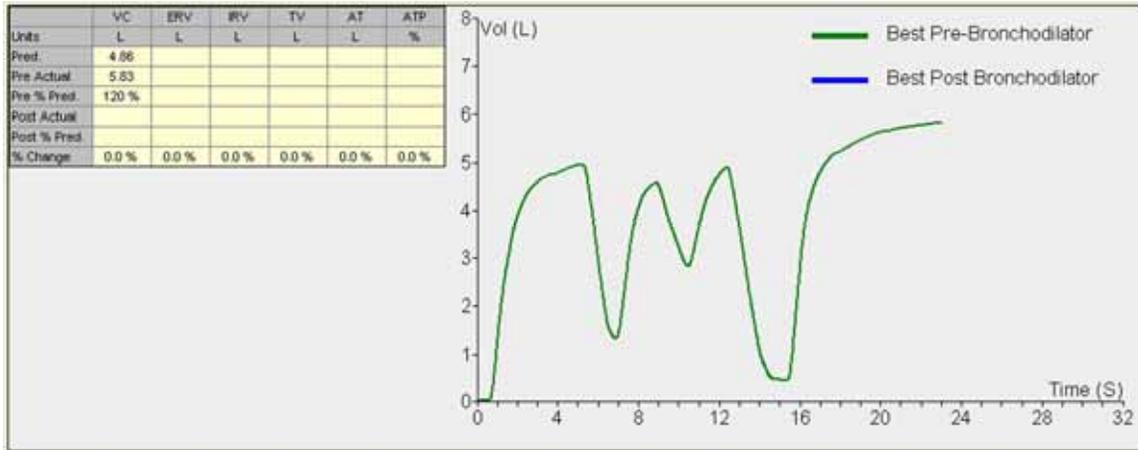


Figure 10
Example of VC test results

Maximum Voluntary Ventilation Test

The maximum voluntary ventilation (MVV) test requires the patient to take several full, deep breaths as quickly as possible for about 12 seconds. This is a demanding test and can lead to hyperventilation. An example of MVV test results is shown in Figure 11.

The MVV test is not normally required in the clinical situations where the AFHCAN Cart is usually installed. Additional details on performing this test are included in the Midmark Operation / Service Manual.

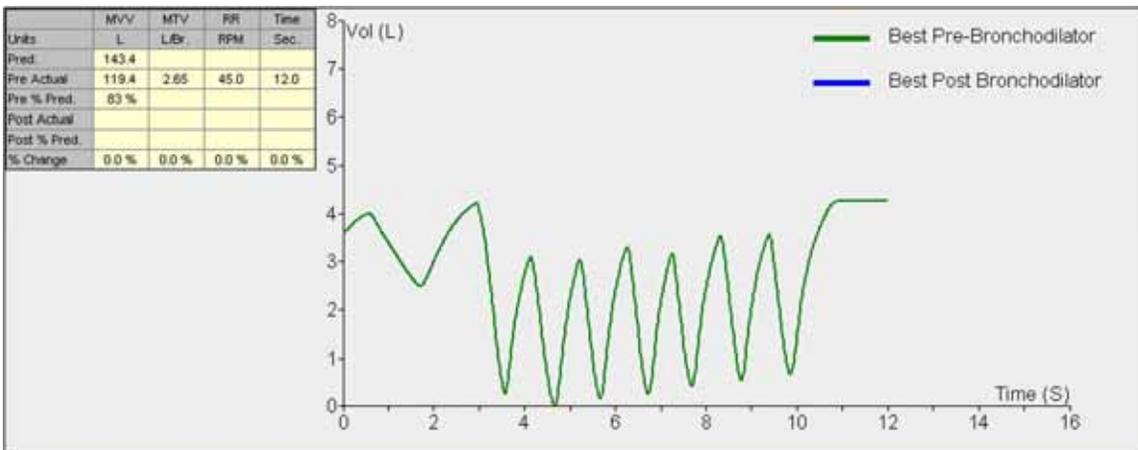


Figure 11
Example of MVV test results

1.1.5.2 Obtaining Spirometry Data

Once all the required patient data has been added, pressing **Start New Report** (or **Continue Report**) from the *Spirometer* screen brings up the Midmark data acquisition screen (Figure 12).

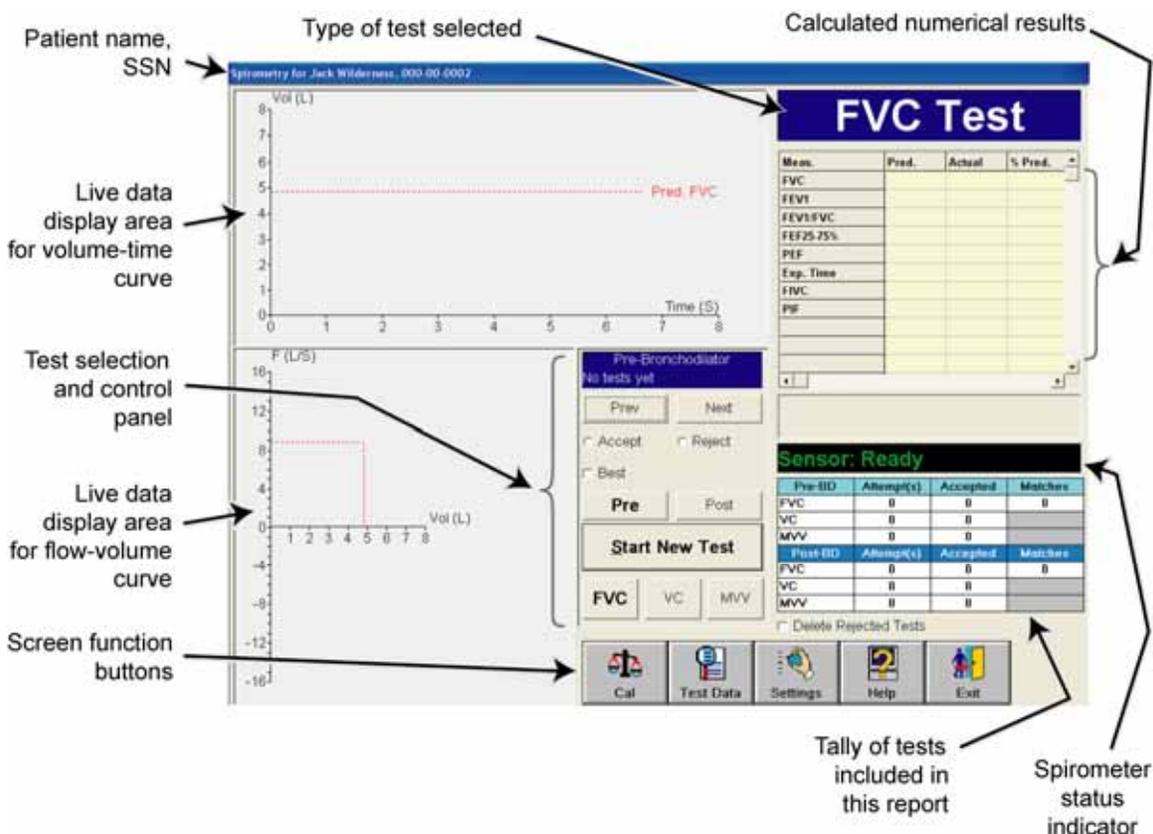
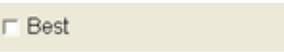


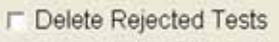
Figure 12
Midmark data acquisition screen

Table 2 describes the items included in the test selection and control panel. Table 3 provides a comprehensive list of the screen function buttons for all Midmark screens encountered in the software.

Table 2
Items in the test selection and control panel

Item	Description
	These three buttons select the type of test. FVC is the default and is the most common spirometry test. The selected test is shown with a bold-face button label. The selected test is also indicated in the upper right corner of the data acquisition screen.

Item	Description
	<p>Before a test is undertaken, these buttons allow the user to identify the test as pre-bronchodilator or post-bronchodilator. The best attempts from each type of test can then be compared side-by-side. The selected condition is shown with a bold-face button label.</p> <p>After some tests have been collected, these buttons select which type of tests can be reviewed in the live data display areas (by using the Prev and Next buttons described below).</p>
	<p>This button begins the selected test. The patient and spirometer must both be ready. When the button is pushed, a message will appear saying Zeroing sensor, please wait. If the sensor is moving or there are any air currents in the vicinity, a message may appear saying Cannot zero spirometry sensor. Hold spirometer still and away from any air currents (usually just to the side of the patient's mouth). Once the sensor has been zeroed, a screen displaying a patient incentive graphic will come up, whereupon the patient can begin the test.</p>
	<p>After some tests have been collected, these buttons sequence through the tests to display them in the live data display areas.</p>
	<p>The information area at the top of the test selection and control panel indicates which type of test is underway or is being viewed. The second line indicates which test in a sequence of tests is underway or being viewed, and the total number of tests of the indicated type that have been taken thus far.</p>
	<p>These selectable items allow the user to choose whether a given test is to be accepted or rejected. The software evaluates the quality of a patient's effort, and will display a message indicating either that the effort was a good one or advising how it can be improved. Any test that rates less than Good test! will be rejected automatically by the software. In the review of one or more test results, the user has the option to manually select whether a given test will be accepted or rejected.</p>
	<p>The software automatically selects one of the tests as the best available effort. In reviewing a sequence of tests, the user may overrule the software's choice and designate a different trial as the patient's best effort. It is the best efforts of pre- and post-bronchodilator trials that are compared side-by-side.</p>


Item	Description
	<p>Saving a set of tests creates a spirometry report. A thumbnail of the spirometry report is displayed on the AFHCAN Spirometer screen. Clicking on the thumbnail brings up the Midmark data display screen (described later). Only accepted tests will be available for display on the data display screen.</p> <p>However, from the data display screen, it is possible to reopen the data acquisition screen. If Delete Rejected Tests has not been selected, any previously rejected reports are still available and can be displayed on the data acquisition screen. The test's status can also be changed to Accepted, making it available for display on the data display screen.</p> <p>On the other hand, if Delete Rejected Tests has been selected (as indicated by a checkmark), rejected tests will be deleted when the report is saved, and cannot be recovered.</p> <p>Once a case has been sent, it is no longer possible to open the data acquisition screen from the data display screen, which means rejected tests can no longer be viewed or recovered regardless of whether this item was selected or not.</p>

Note: The screen function buttons in Table 3 may be disabled or replaced with other buttons from time to time depending on what the program is doing. When a button's function is not available, the button is grayed out.

Table 3
Midmark screen function buttons and applicable screens

Button	Screens ¹	Function
	A	Brings up the Midmark calibration screen.
	A	Brings up a dialog box where the medical indication, bronchodilator used, barometric pressure, ambient temperature, and dry air conditions can be entered or selected. The bronchodilator can also be selected from the data display screen. The parameters entered will appear on a report printout, if one is generated.
	A C	Pressing the Settings button opens up a nine-tab dialog box. Details on settings are provided in Section 1.1.5.5. The data display screen can access the settings dialog box with the alternative form of the settings button shown at the right. Parameters set in the Settings dialog box persist in the memory and become the default settings for future sessions.

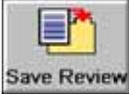

Button	Screens ¹	Function
	A D C	Pressing the Help button brings up the Midmark help pages. Not all features described are available in the AFHCAN installation.
	A D C	Pressing the Exit button exits from the current screen and returns to a previous screen in the sequence.
	A	When one or more tests have been performed, this button replaces the Exit button on the data acquisition screen. Pressing this button saves the report and returns to the Spirometer screen.
	D	The File button appears when viewing an existing report. Pressing the File button brings up a grayed-out menu for archiving or emailing the report. These functions are handled by the AFHCAN software.
	D C	The Print button prints a spirometry report. The print format can be modified in the Settings popup window (see Settings below).
	C	This button allows a quick check of the calibration to be done. It only requires one pump of the calibration syringe. (If you are uncertain whether this is adequate for your purposes, perform a full calibration.)
	C	Initiates the sequence for performing a full calibration, and activates the Start Cal button.
	C	This button brings up the calibration pump window, which prompts the user to pump the calibration syringe, and provides a pump rate guidance graphic, which helps the user gauge the rate at which the syringe is pumped.

Note 1: the *Screens* column indicates which screens the button can appear on, as follows:
A = data acquisition screen; D = data display screen; C = calibration screen.

Patient Incentive Displays

When the **Start New Test** button is pushed on the data acquisition screen, the Midmark software brings up a patient incentive display, which is designed to help patients make their best efforts. There are five patient incentive displays that can be selected via the Settings Dialog Box, Incentives/Miscellaneous tab (see 1.1.5.5 for information). Patient incentive displays respond to the patient's level of effort as if the patient is blowing on the display directly (blowing out candles, for example). Three of the five available displays are shown in Figure 13.

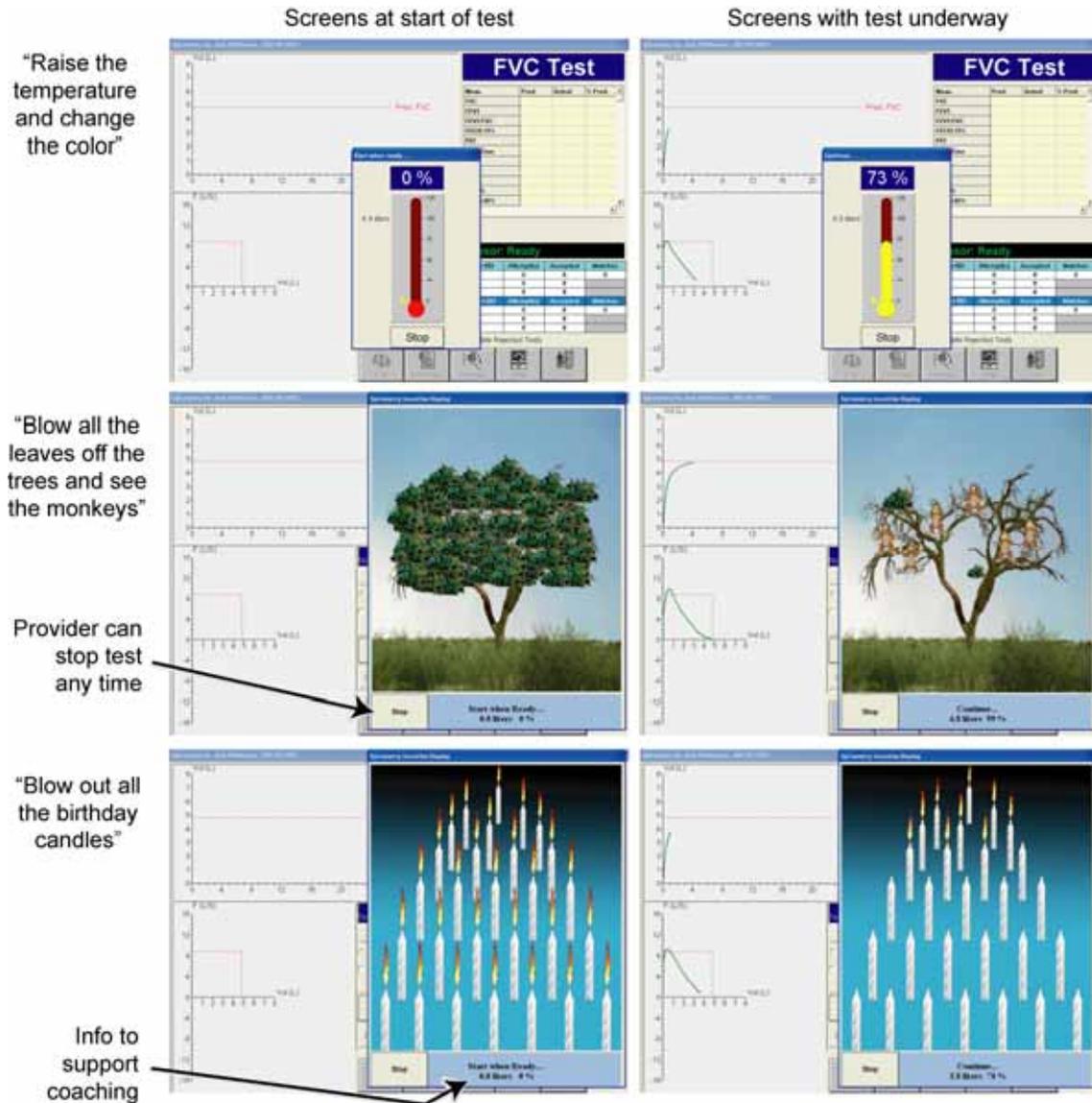


Figure 13
Patient incentive displays

Patient Feedback Messages

If the patient inhales through the mouthpiece, the inhalation will be recorded and the test will end automatically. If all that is required is the FVC test, then the operator can stop the test when it is obvious the patient has emptied his or her lungs. To stop the test, press the **Stop** button on the patient incentive display.

On completion of the test, a message providing feedback on the patient's performance will be displayed. The goal is to obtain the "Good test!" message shown in Figure 14. Also shown is an example offering advice on how to improve performance. Details on all the possible messages can be found in the Midmark Operation / Service Manual.

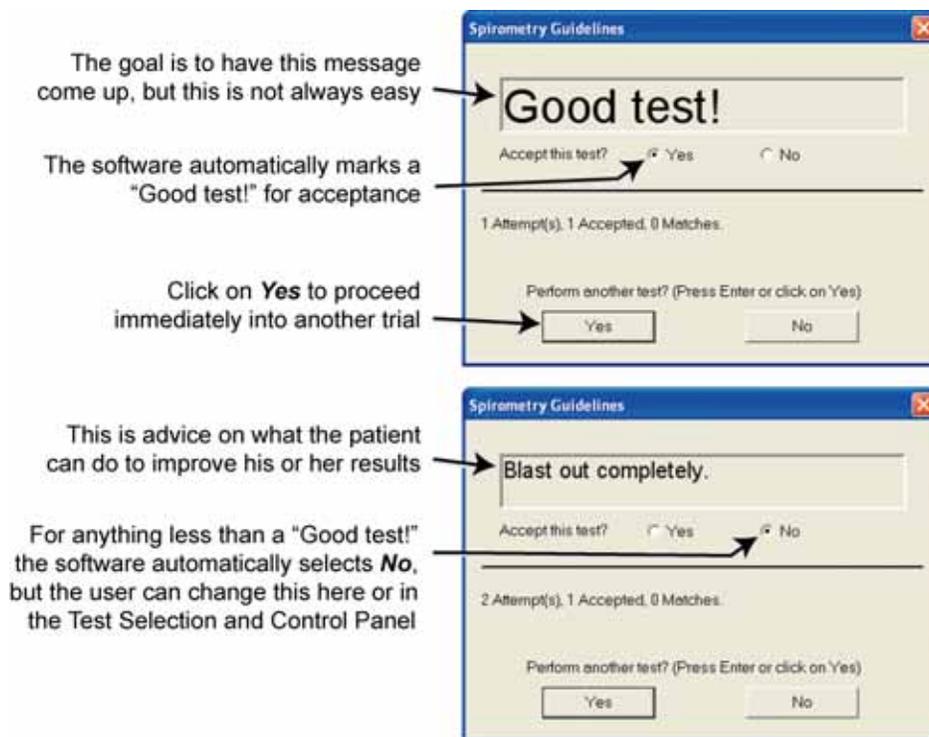


Figure 14
Examples of patient feedback messages

When all desired trials and tests have been obtained, click on **No** in the feedback window, then **Save Review** in the data acquisition screen. A thumbnail of the spirometry report will now appear in the left bar area of the *Spirometer* screen.

1.1.5.3 Viewing Spirometry Data

Clicking on a spirometry thumbnail in the *Spirometer* screen (or in the *Case* screen) brings up the Midmark data display screen. Initially the screen display is in the data format, which is the first tab selection (*Summary*), as shown in Figure 15.

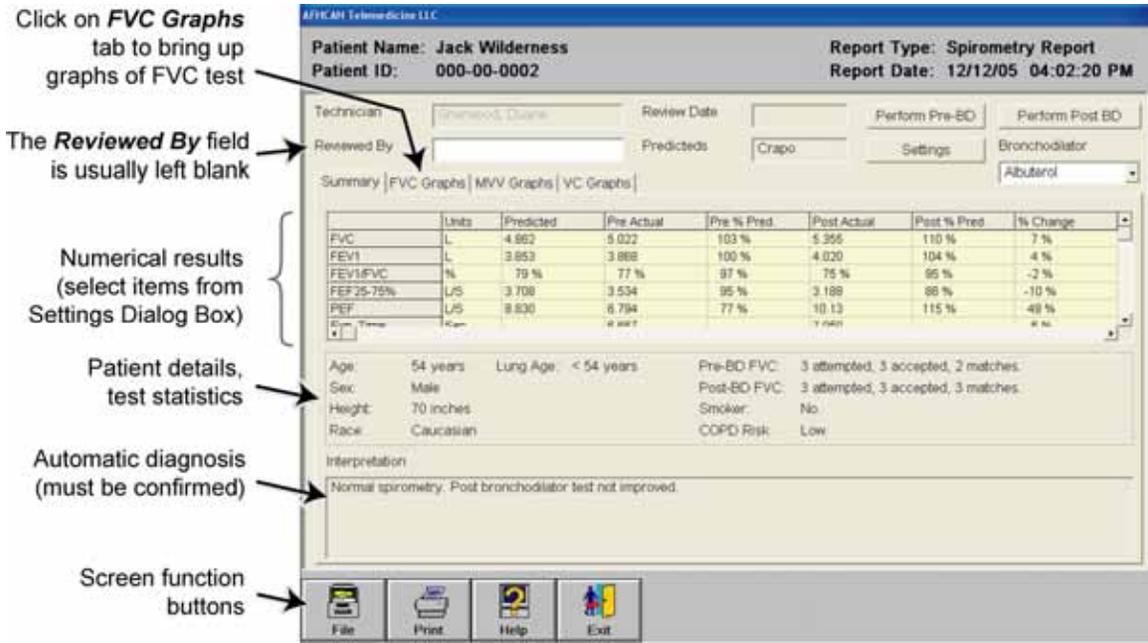


Figure 15
Midmark data display screen – numerical data

Clicking on the **FVC Graphs** tab brings up a graphical display of the spirometry results (Figure 16).

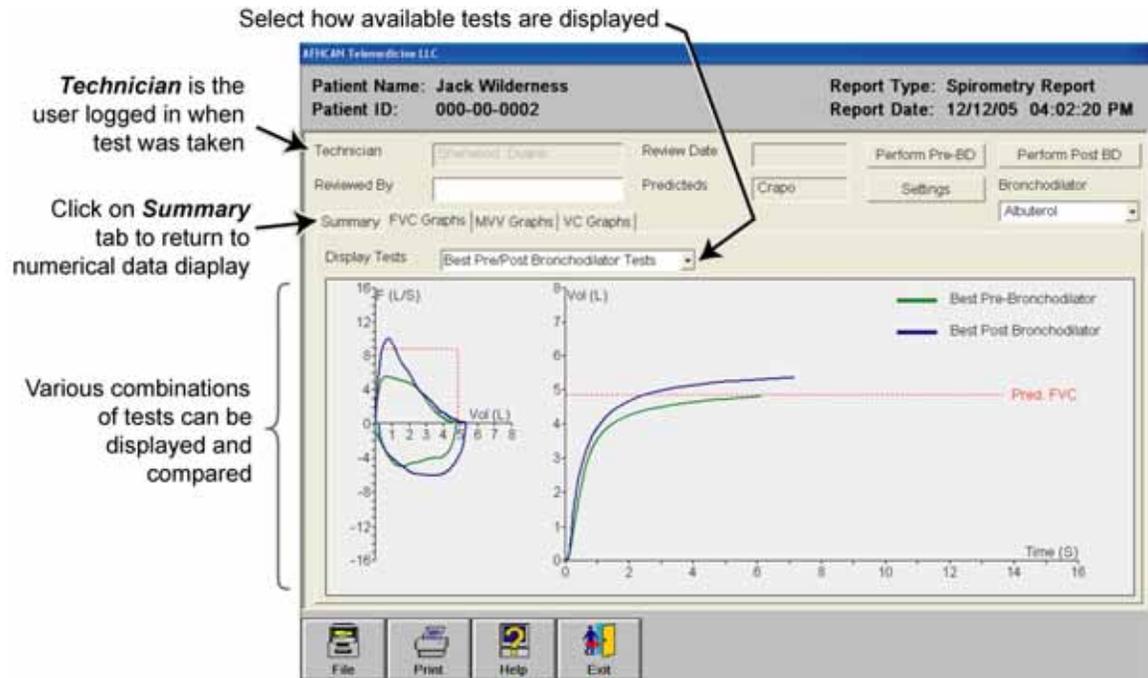


Figure 16
Midmark data display screen – graphical data

Viewing Active Data and Viewing Fixed Data

The Midmark data display screen can be called up in two different situations:

- during the process of adding a new spirometry report in a case (data is active)
- when viewing an existing report after a case has been sent (data is fixed)

When a new spirometry report is being added to a case, it opens up what can be called a "spirometry session." The session remains open until the case is closed and sent (this includes placing the case on hold). At the beginning of the session, supplemental patient information must be added before any spirometry tests can be performed. As long as the session is open, the following points apply:


- the supplemental patient information entered at the beginning remains available to the Midmark spirometry software
- one or more separate spirometry reports can be created
- one or more trials can be added to a given report


These points hold true even after a spirometry report has been saved to a case. If the spirometry session is still open (in other words, new spirometry data has been added, but the case has not been closed and sent), it is possible to open a spirometry report from the **Case** screen and add data to it.

You can tell whether a spirometry report is still active or not by looking at the Midmark data display screen. Table 4 describes the features of the Midmark data display screen which are available when a spirometry session is open, but which are unavailable after the data has been sent.

Table 4

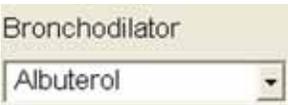Data display screen items that are active when a "spirometry session" is open

Item	Description
	Clicking on this button brings up the Midmark data acquisition screen with Pre already selected. New pre-bronchodilator spirometry tests can be added to the report.
	Clicking on this button brings up the Midmark data acquisition screen with Post already selected. New post-bronchodilator tests can be added to the report.
	When this area is active, it is possible to identify or change the bronchodilator used in the test.
	Brings up the nine-tab settings dialog box. Functional equivalent of the Settings function button. 

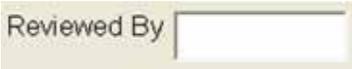Item	Description
	<p>Normally this field is left blank by the originating provider. In the Midmark software, the Reviewed By field would be used by the healthcare provider responsible for reviewing the case and affirming or modifying the diagnosis. Once the case is sent, however, this field is locked. The reviewing provider would normally enter a diagnosis or instructions in a case comment, not in the spirometry report itself.</p>
	<p>This field is not directly accessible by the user. The Review Date field works in conjunction with the Reviewed By field. If for some reason there is information entered in the Reviewed By field, then the date the information was entered would be displayed here by the Midmark software.</p>


1.1.5.4 Calibration Screen

With the Midmark Digital Spirometer, calibration confirms that the instrument is operating properly. Calibration is done using the Midmark 3-liter calibration syringe. A mouthpiece (preferably one reserved for calibration) is placed on the spirometer handle, and then it is fitted into a rubber adapter on the syringe. Following directions on the screen, the user pumps the syringe a number of times until a calibration result is displayed. (Use a helper to hold the spirometer in place.) If the result is accepted, the calibration data is stored in a thumbnail and is available for inclusion in a case. The calibration results will be stored and used by the Midmark software until the next time calibration is performed.

The American Thoracic Society recommends that calibration be performed once each day (assuming of course that the device is being utilized that day). It is a good practice to maintain records of when calibration was performed. One option is to create a fictitious patient named "Spirometry Calibration," and to save (and archive) all calibrations to test cases associated with that patient name. That way there is a permanent record of calibrations that anyone can easily access. Alternatively, a calibration can be saved into a real case along with the patient's spirometry data, or your facility can print and file calibration reports.

The Midmark software allows the user to specify the calibration expiration interval. Using the Settings dialog box, intervals from 1 to 30 days can be entered. Each time spirometry is done, the software will compare the date of the last calibration against the calibration expiration interval. If the calibration has expired, the user will be notified and invited to perform a calibration before proceeding.

The most common way to initiate calibration is by pressing the **Calibrate** button on the AFHCAN *Spirometer* screen. Doing so brings up the screen shown in Figure 17.

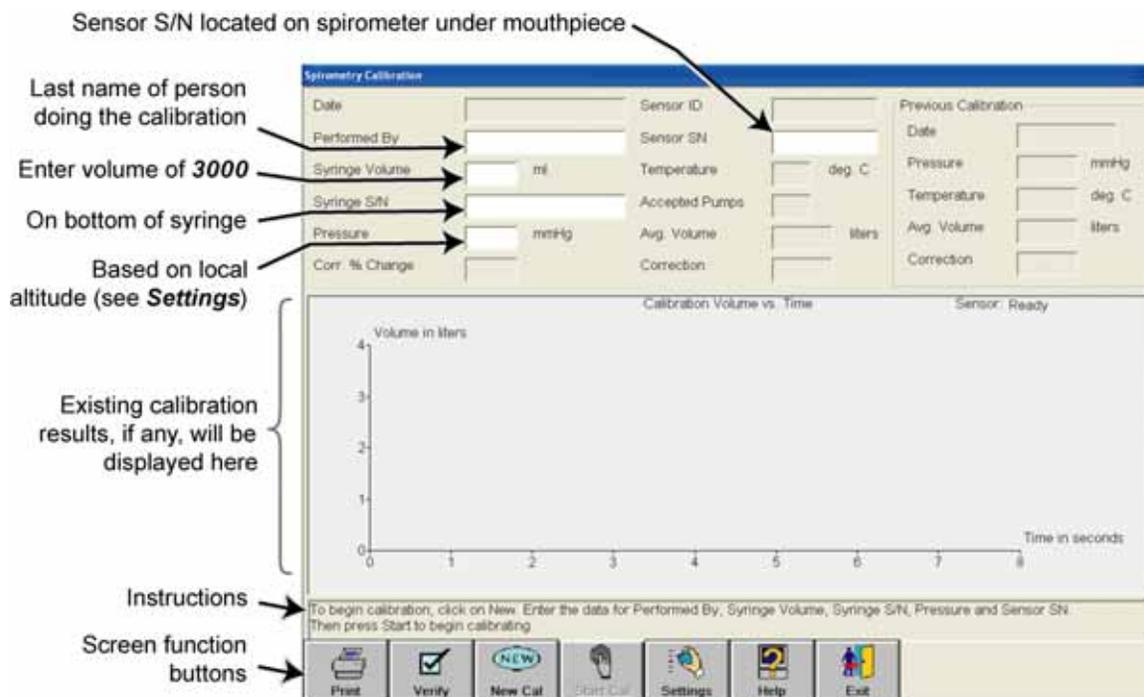**Figure 17**

Midmark calibration screen – first calibration

Figure 17 shows how the screen appears the very first time a spirometry report is created at an AFHCAN Cart. The blank fields must be filled in. The callouts on the illustration show where to obtain the desired information.

The **Pressure** entry requires a little more explanation. The following are options for entering or obtaining an appropriate value:

- if you happen to have a reading of local atmospheric pressure, say from a daily weather report, and the units of pressure are in millimeters of mercury (mm Hg), then you can enter that number directly
- if the elevation of your location is within a few hundred feet of sea level, you can enter the value for standard atmospheric pressure of 760 mm Hg
- if you know the elevation of your location (in feet above sea level), then you can allow the Midmark software to computer a standard reference pressure for that elevation as follows:
 - go to the **Settings** dialog box
 - select the **Calibration** tab
 - select **Calculate Barometric Pressure from Altitude**
 - enter the elevation of your facility in the **Altitude Above Sea Level** text box

When the Midmark software calculates a reference pressure, that revised pressure will appear in the Pressure text box in the calibration screen.

The accuracy of the Midmark Digital Spirometer is not affected significantly by daily fluctuations in atmospheric pressure. However, if the device is moved to a new location, and the difference in

elevation is significant (say, 1,000 feet or more), then a new reference pressure should be obtained, and a new calibration should be performed.

Once the instrument has been calibrated one time, the Midmark calibration screen will come up with the fields already filled in, as shown in Figure 18.

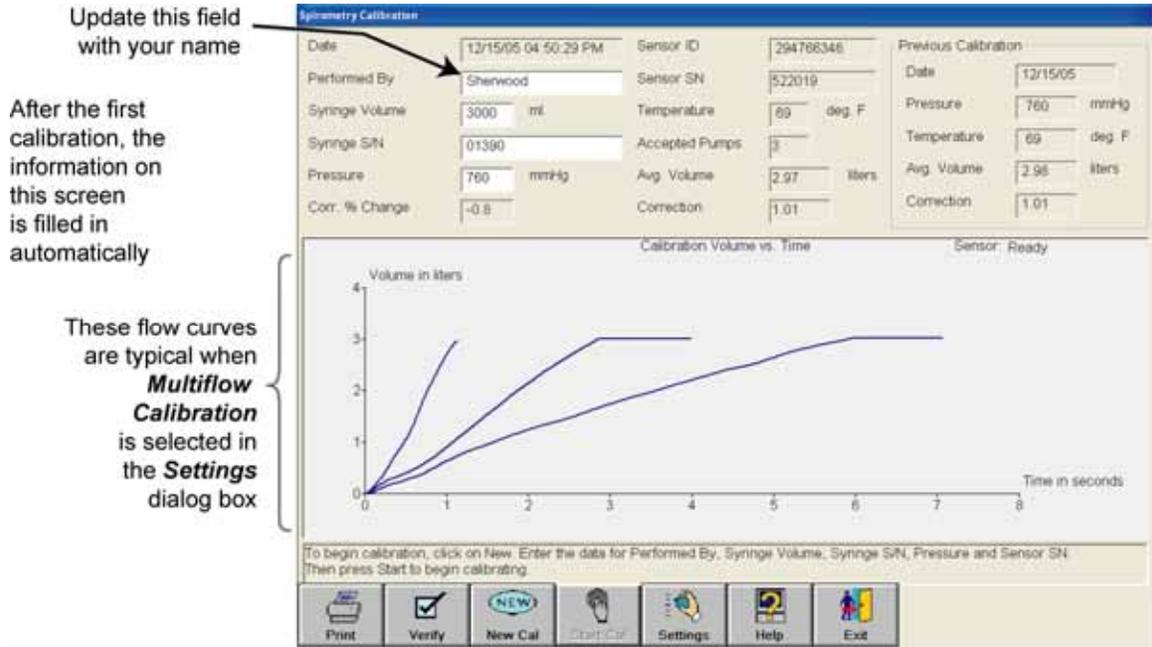


Figure 18
Midmark calibration screen – subsequent calibrations

Calibration Pump Window

Pressing **New Cal** in the calibration screen begins the process of creating a new calibration report. At this point you can update the various fields on the screen.

Note: The name displayed in the **Performed By** field comes from the previous calibration report, not from the name of the current **User** (the person logged into the AFHCAN software). The person actually doing the calibration should enter his or her last name in the **Performed By** field.

When the data is satisfactory, place an unused mouthpiece on the spirometer, place the adapter on the calibration syringe, ensure the plunger has been pushed all the way in, and press the **Start Cal** button. Pressing the **Start Cal** button will bring up the calibration pump window (Figure19).

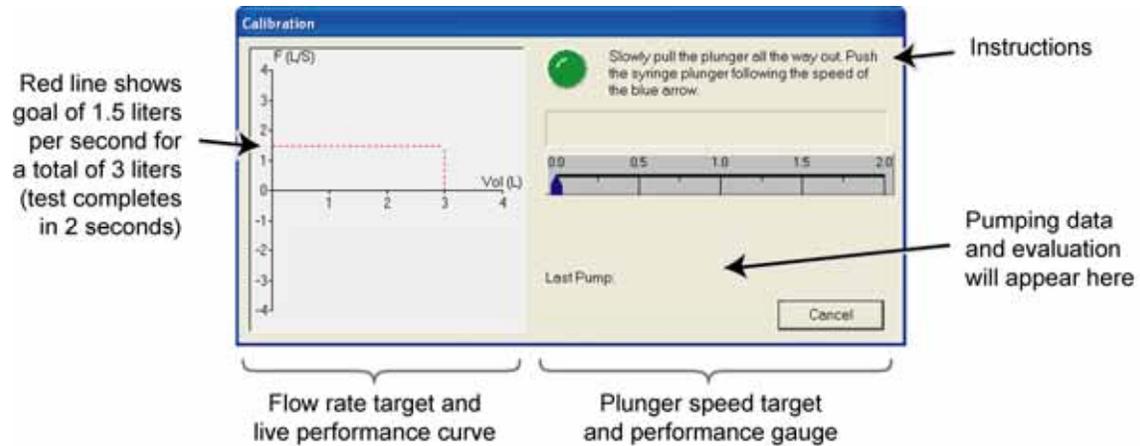


Figure 19
Calibration pump window – single flow

Follow the instructions on the calibration pump window and slowly draw the pump handle out all the way. As soon as the pump handle begins to push back in, the display will go live. On the left side, a trace of the flow rate will appear. On the right side, the blue pointer will start moving and a performance bar will appear. The goal is to have the performance bar stay even with the moving pointer. It takes a little practice to get used to it. There is no harm in doing several calibrations in a row. Only the most recent one is saved in the Midmark software. Practice calibration reports can be unchecked so they are not saved in a case.

Note: It is important to draw the plunger out slowly. If it is drawn out too fast, it will bounce a little bit when it reaches its limit of travel. The spirometer is sensitive enough to detect the bounce as forward flow, whereupon it begins moving the needle to the right. If you are not actually ready to begin pushing yet, it can be hard to catch up so the performance bar stays even with the pointer. If needed, simply restart the test.

Figure 19 shows the calibration window that comes up when Multiflow calibration testing has not been selected. Figure 20 shows how the calibration pump window looks when Multiflow calibration testing has been selected, with a calibration test underway.

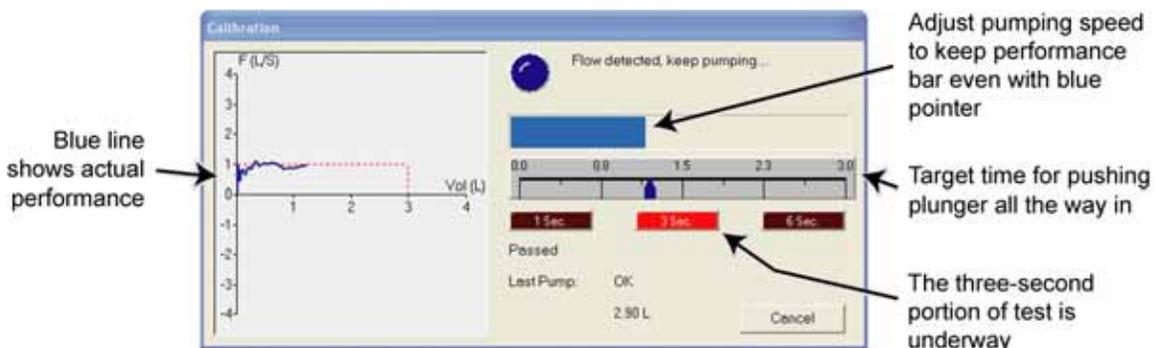


Figure 20
Calibration pump window – multiflow test underway

1.1.5.5 Settings Dialog Box

The Settings dialog box allows the user to control the spirometer's performance, alter the default display and print formats, and select system preferences. Figure 21 shows the basic plan of the window. Click on tabs at the top of the window to select an individual tab page. If any settings are changed, click on Apply at the bottom of the window to activate the setting. Click on **OK** to save all settings and exit. Selections will be saved in memory and will become the default settings for subsequent spirometry reports. Click on **Cancel** to cancel any changes not yet applied and exit.

Figures 22 through 24 show more detailed images of selected tab pages as well as typical settings for items on that page. For more information on items not covered in this manual, refer to the Midmark Operation / Service Manual.


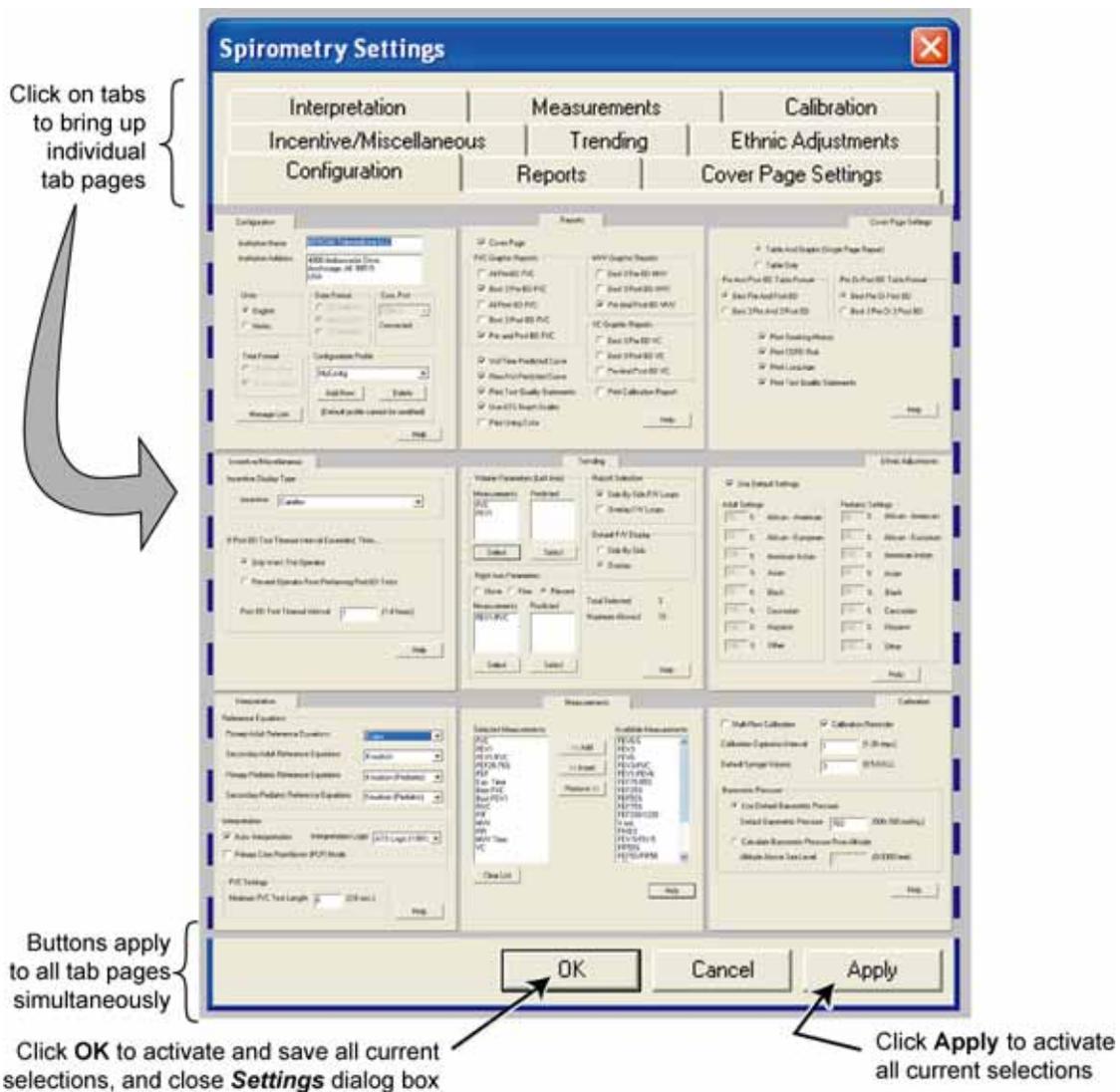


Figure 21
Available tabs in settings dialog box

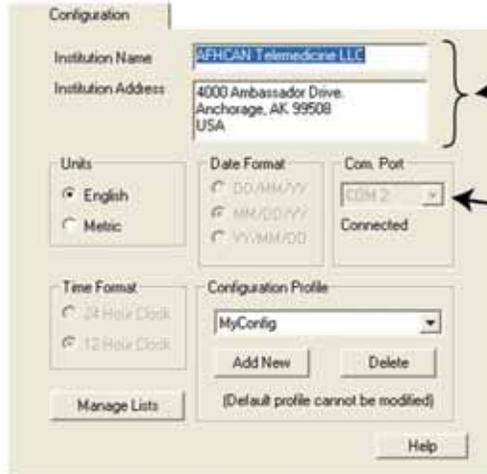
The **Configuration** tab is for some general configuration items



Information entered here will appear at the tops of printed reports

A serial COM port is not applicable to USB-type spirometers

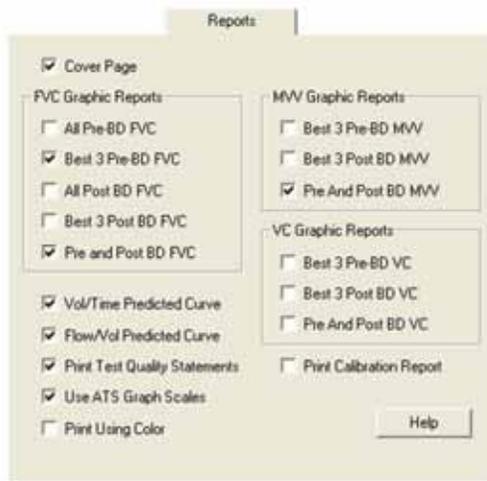The **Reports** tab allows you to select items to appear on printed reports




The **Cover Page** tab allows you to select items to appear on the cover page of printed reports



Figure 22
Settings tabs – first row

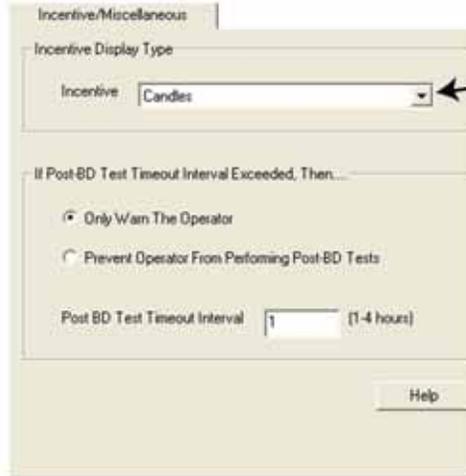The **Incentive/Miscellaneous** tab allows selection of patient incentive displays and Post-BD warning and timeout




Select the patient incentive display that works best for you and your patient

Since post-BD tests are usually done a short time after administration, these user convenience features are of limited value

The **Trending** tab applies to trending reports, which is not available in the AFHCAN installation



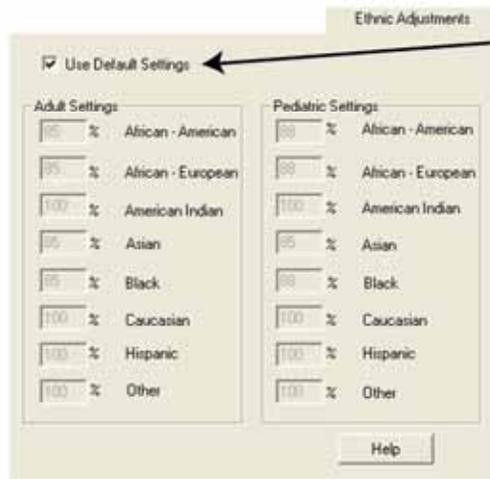
The **Ethnic Adjustments** tab allows factors used in calculating predicted performance to be changed for various ethnic groups



It is best to simply leave this checked


Figure 23
Settings tabs – second row

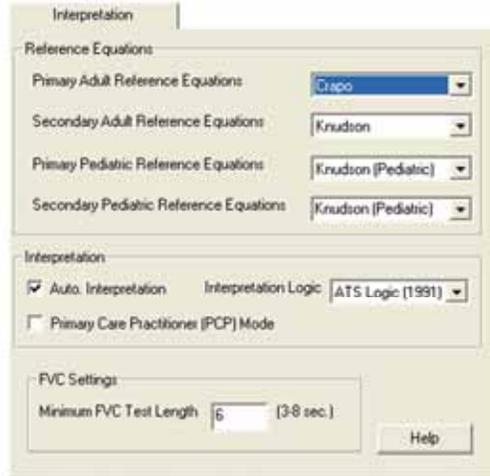The **Interpretation** tab allows selection of reference equations used in predicting performance and interpreting results



The **Measurements** tab allows various calculated values to be included on the report

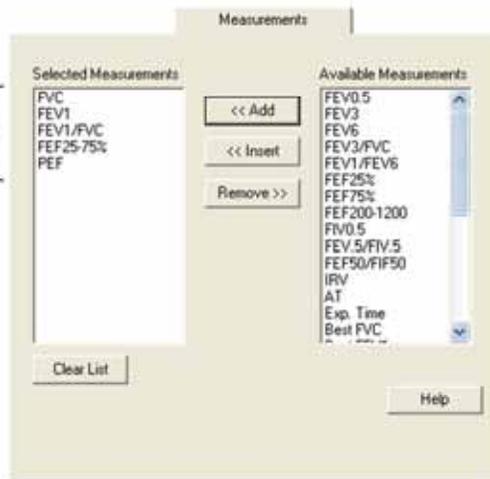Standard items to include on report




It is possible to add or delete items based on your facility's procedures or the consulting provider's request (default values are recommended for most situations)

The **Calibration** tab allows various calibration settings to be controlled

Multi-flow calibration works with Calibration Pump Window (required if used for OSHA testing)

If local barometric pressure is unknown, click here

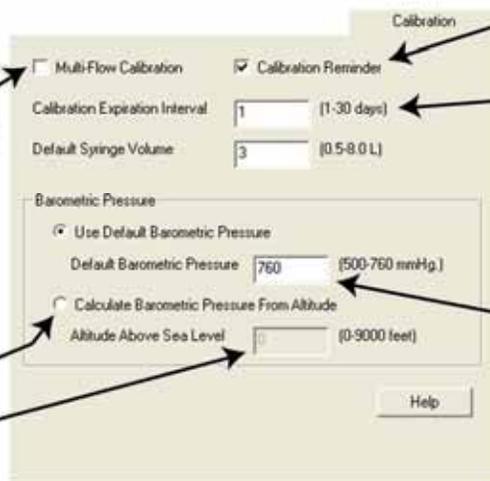.... and enter the approximate elevation of your facility here



When selected, reminds user that calibration is due

Recommended frequency for performing calibrations (the ATC recommends daily calibrations; your facility may have a different requirement)

If known, enter a typical barometric pressure from a daily weather report

Figure 24
Settings tabs – third row

1.2 Particulars of the AFHCAN Installation

The AFHCAN software serves as a host program for the Midmark software. This means the AFHCAN software calls up the Midmark software on an as-needed basis, and works with the Midmark software to provide full functionality. Functions such as filing and archiving reports, for example, are handled by the AFHCAN software, so the corresponding features in the Midmark software have been disabled. Because of technical requirements in the software, various other features described in the Midmark Operation / Service Manual are not available in the AFHCAN installation. None of the disabled features compromise the intended functionality. Buttons, text entry boxes, or other items are typically grayed out if they are not available.

The AFHCAN software provides some information that is required by the Midmark software in order to work (e.g., patient, gender, age,). Integration with the AFHCAN Cart and software means there is no setup or installation required.

1.2.1 Version Compatibility

When sending a case, spirometry reports are not backwards compatible with older versions of the AFHCAN software. This means a report created using AFHCAN Client software Version 4 cannot be sent to providers using the earlier Version 3.4. The AFHCAN software will notify the sender if this situation occurs, and will prevent the case from being sent. It is always possible to print out the spirometry report, scan it into a new case, and send that case.

1.3 Warnings and Cautions

Refer to the Midmark Operation / Service Manual provided with the device for warnings and cautions.

In additional, please observe the following points:

- any automatic analyses provided by software must be confirmed or revised by a qualified physician
- do not use the spirometer with patients who are nauseated or in situations where liquids can enter the mouthpiece and flow into the pressure ports.
- be aware of the patient and discontinue the test if any abnormality is observed (such as hyperventilation)
- discard the disposable mouthpiece after each use
- use a separate mouthpiece for calibration

Section 2 – Operation

2.1 Basic Operating Procedures

Refer to Section 1 for descriptions and illustrations of the controls, indicators, screens, and other items referenced in these procedures.

2.1.1 Starting the Spirometer

Proceed as follows:

1. Plug the USB connector into an available USB port on the Cart or computer.
2. From the **Add To Case** screen, press the **Patient** button to identify the patient (search, add, or edit patient information as necessary).
3. From the **Add To Case** screen, press the **Spirometer** button to begin a spirometry session. This brings up the **Spirometer** screen.

2.1.2 Calibration

With the Midmark Digital Spirometer, calibration is more of a functional checkout than an actual adjustment to the instrument. The frequency with which you perform calibration depends on your facility's policies and requirements. The American Thoracic Society recommends daily calibration. A reminder interval for doing calibrations can be entered under the **Calibration** tab in the Settings dialog box.

To maintain a consistent and readily accessible record of calibrations, consider creating a test case with a patient named "Calibration Reports" or something similar. All cases for this patient would then contain the calibration reports. The cases can be archived and then can be easily searched. As an alternative, you can include the calibration report along with the spirometry reports for a real patient.

Proceed as follows:

1. From the **Spirometer** screen, press the **Calibrate** button. This brings up the **Spirometry Calibration** screen. After a moment, the previous calibration results will be displayed, and the sensor status will be indicated at the right of the screen.

Note: The **Spirometry Calibration** screen can also be accessed from the data acquisition screen by pressing the **Cal** button.

2. Press the **New Cal** button.
3. Confirm or enter the necessary data in the data entry fields at the top of the screen (if you are doing the calibration, enter your last name in the **Performed By** field).
4. Press the **Start Cal** button. The **Calibration** pump window will come up.

5. Have someone assist by holding the spirometer in the adapter while you pump the plunger, as shown in Figure 25. (This is an excellent way to involve the patient in the process.)



Figure 25

Double-teaming the calibration syringe

6. Follow the instructions displayed in the window. As soon as the spirometer detects forward flow (plunger moving in), it will start moving both the pointer and the performance bar. Adjust pumping speed as needed to keep the performance bar even with the pointer.

Note: When pulling the plunger out, do so slowly and gently. If pulled out too quickly, the plunger may bounce a little on reaching the limit of travel. This bounce may be enough to trigger the performance bar. If you are not ready, it may be hard to catch up to the pointer. At worst, this may mean the system asks you to repeat the attempt.

7. Pump as many times as requested, each time matching the speed of the pointer. Follow any instructions that may appear. The system will reject any unsatisfactory attempts, inform you of why, and will ask you to repeat the trial.
8. After a sufficient number of satisfactory attempts have been made (usually 4 or 5), the screen will display the *Verification Pump* window. If the indicated results are within the ranges listed, click **Yes** to accept the calibration.
9. Click **Exit** to return to the *Spirometer* screen. Each calibration report generated will appear as a thumbnail in the left bar area. A green checkmark will appear automatically indicating that these reports will be saved to the case.
10. To uncheck any reports that you wish to discard, click on the green checkmark next to those reports. The green checkmark will disappear and that report will not be saved to the case. Click the box again to bring the checkmark back.
11. If this case has been set up to record only calibration results, click **Save** to save the reports and return to the *Case* screen; process the case as desired (e.g., send to archives). Otherwise, proceed with obtaining a spirometry report.

2.1.3 Obtaining a Spirometry Report

2.1.3.1 Patient Preparation and Coaching

The quality of the spirometry report is directly related to the patient's level of effort. Being clear about how to perform the test, and conveying enthusiasm to the patient, are very helpful in ensuring the best effort is obtained.

Proceed as follows:

1. Explain that the goal is to inhale deeply and to exhale forcefully, quickly, and fully.
2. Using your own clean mouthpiece, demonstrate how to perform the test.
3. Explain that teeth and lips go outside the mouthpiece.
4. Give the patient a clean mouthpiece, and let him or her blow through it (handle not attached).
5. Loosen or remove any tight clothing.
6. Remove any foreign objects or dentures from the patient's mouth.
7. Let the patient know he or she can stop the test if feeling dizzy, ill, or uncomfortable.
8. If pre- and post-bronchodilator tests will be performed, have the medication handy.

2.1.3.2 Conducting the Test

Proceed as follows:

1. On the *Spirometer* screen, click on the **Start New Report** button.
2. If this is the first test of the session, enter the required fields for supplemental patient data (race, height, weight, smoking history, and medications), and click **Continue Report**.
3. From the Midmark data acquisition screen, click **Test Data** to enter information about symptoms and the bronchodilator used, if desired.
4. Ensure the type of test selected is **FVC**.

Note: Because **VC** and **MVV** are not normally done in the settings where the AFHCAN Cart is deployed, the **FVC** test is assumed.
5. If this is a pre-bronchodilator test, ensure that **Pre** is selected (the selected button is shown in bold letters). If you have already performed the pre-bronchodilator tests, and this is a post-bronchodilator test, select **Post**, and enter the bronchodilator used in the pull-down list in the upper right area of the *Spirometer* screen.
6. Place the patient's mouthpiece on the handle ensuring it is fully seated in the slot.
7. Place a disposable nose clip on the patient's nose (highly recommended).
8. Have the patient hold the mouthpiece still beside his or her cheek.
9. Press the **Start New Test** button. An incentive display will come up.

10. Tell the patient:
 - a. inhale as deeply as possible
 - b. place mouthpiece in mouth with teeth and lips over outside and lips making a good seal
 - c. blow out as fast, hard, and long as possible, at least 6 seconds (but not more than 15 seconds)
 - d. inhale deeply (optional for full flow volume loop)
11. If step *10-d* is not required or not done, press **Stop** on the incentive display to end the test.
12. A *Spirometry Guidelines* window will come up providing feedback on the test and guidance on what the patient can do to obtain a **Good Test!** message.
13. If anything but the **Good Test!** message appears:
 - a. the test will automatically be marked **No** for accepting the test; you have the option to select **Yes**, which means the results for this attempt will be saved in the data included in this report
 - b. use the information displayed to coach the patient on how to improve performance

Note: The **Good Test!** message is fairly difficult to obtain. Let the patient know this, and encourage them to keep up the good effort.
14. The *Spirometry Guidelines* window allows you to select whether to perform another test of the same type or not:

Note: A minimum of three tests should be done, but no more than eight tests.

 - a. if you select **Yes**, repeat steps 9 - 14
 - b. if you select **No**, the Midmark software returns to the data acquisition screen, which still allows you to initiate further testing
15. On completion of the testing, click on the **Save Review** button to close the Midmark software and return to the AFHCAN *Spirometer* screen.

2.1.3.3 Reviewing and Selecting Tests

The Midmark data acquisition screen allows you to review a series of tests, decide which to keep, and to manually identify the patient's best effort. This is done using the test selection and control panel in the middle of the screen. Proceed as follows:

1. To review the newly generated tests, click on **Prev** or **Next**. The software will display each test of the given type in sequence.
2. To switch the type of test being reviewed, select **Pre** or **Post** (for pre-bronchodilator tests and post-bronchodilator tests).
3. As each test is displayed, you can select whether to **Accept** or **Reject** the test, or you can indicate that the test being displayed was **Best** (the patient's best effort).

Note: Rejected tests can be recovered as long as the current session remains open, but are lost once the case itself has been saved.

4. To take additional tests, click on **Start New Test**.
5. On completion of the testing and review, click on the **Save Review** button to close the Midmark software and return to the *Spirometer* screen. The spirometry report just completed will now appear as a thumbnail in the left bar area. The green checkmark will automatically appear in the box to the right, indicating that this report will be saved to the case.

2.1.3.5 Saving the Spirometry Report to a Case

From the *Spirometer* screen, click the **Save** button. This returns to the *Case* screen. At this point, the data has been saved to the case, but if the report is opened again before the case is sent, data from more tests can be added.

It is a good practice to describe the quality of the patient's efforts. Were they 100% cooperative and making their best effort? Were there distractions that could have affected performance. On the *Case* screen, press the **Add Comment** button to enter this kind of descriptive information into the case. If possible, try to correlate the descriptions with individual tests.

Once the case is sent, the data in the spirometry report becomes part of a permanent medical record and cannot be changed. New spirometry reports could be added to the case, however.

2.1.4 Viewing Spirometry Reports

Proceed as follows:

1. From the *Spirometer* screen, or from the *Case* screen, click on the thumbnail for a spirometry report. This will open the Midmark data display screen.
2. On first opening the data display screen, the *Summary* tab will be active and numerical results will be displayed. Scroll up or down, and left or right, to bring data into view.
3. To see graphics of the patients performance, click on the *FVC Graphs* tab. (If other types of spirometry tests were done, click on the applicable tab.)
4. To display various combinations of graphs, click on the *Display Tests* pull-down window.
5. If the data is active, you can press the **Perform Pre-BD** button or the **Perform Post-BD** button in the upper right corner of the screen, and add more tests to the report.

Note: If the data is not active (in other words, if the case has been sent and the data has been frozen as part of a permanent medical record), these buttons and other features will be grayed out.

6. To print the spirometry report, click on the **Print** button.

Note: Various items in the *Settings* dialog box give the user some control the format and content of the printed reports.

7. To close the report and return to the previous screen, click **Exit**.

Section 3 – Clinical Considerations

These pointers are to be used as a supplement to your clinical experience.

3.1 Guidelines for Clinical Success

- coaching the patient in how to make a "best effort" is important
- demonstrate yourself how to do the test
- have the patient sit up straight
- it is important for the patient to give a good effort
- always use nose clips to help the patient achieve a good seal
- goal is three good tests
- no more than eight attempts

3.2 Common Mistakes

- lips not making a tight seal
- patient coughs during test
- spirometer not calibrated recently (calibrate once a day – not required for every patient)

3.3 Tips and Tricks

- help the patient relax by describing the first attempt as "practice"
- do not attempt spirometry on a patient who is uncomfortable, acutely ill, or uncooperative
- stop the test if the patient becomes tired or uncomfortable (if a good curve cannot be obtained in 8 attempts, allow the patient to try on another day)
- when calibrating the spirometer, have the patient or another helper pump the calibration syringe for you
- give the patient a chance to empty his or her bladder before the test to reduce the risk of leaking urine if the patient coughs
- most children of age seven (and some younger) can be given a spirometry test

Section 4 – Routine Maintenance

4.1 Care and Cleaning

Since the mouthpiece is the single-use, disposable type, cleaning is limited to a general wipe-down of the handle and cable with cloth dampened with a mild detergent. Do not get any moisture on the pressure ports at the top of the handle.

4.2 Replacement Items

The USB version does not require batteries. The battery compartment may have been glued shut.

Box of 25 Disposable Pneumotach Mouthpieces, Midmark part number 2-100-1200

4.3 Elementary Troubleshooting

Table 5 shows offers some troubleshooting guidelines. Additional troubleshooting information can be found in the Midmark manual, or on their website.

Table 5

Troubleshooting guidelines

Problem	Possible Solutions
Midmark software indicates the sensor is not ready or cannot be found	<ul style="list-style-type: none">• ensure the USB connector is connected to a suitable USB port on the Cart or computer• try exiting the Midmark software and coming back into it from the Spirometer screen
calibration results not within a satisfactory range	<ul style="list-style-type: none">• may indicate a problem with the instrument• in general, the pressure ports on the handle are not designed for user cleaning or maintenance; however, it does no harm to inspect the openings to be sure no foreign objects or material is present
cannot zero the spirometer	<ul style="list-style-type: none">• the spirometer is very sensitive; if any air movement is detected, the spirometer cannot zero• while zeroing is underway, keep instrument still and away from any miscellaneous air currents in the environment
other problems	<ul style="list-style-type: none">• check the Midmark Operation / Service Manual• contact AFHCAN Customer Support at 888-449-4435• check the Midmark online help pages by clicking Help on their screens

Index

- Accept radio button
 - function and use, 12
 - on data acquisition screen (image), 11
- adapter
 - for calibration syringe (image), 3
 - general reference, 19, 21
 - in position (image), 30
- Add Comment button, 33
- Add To Case screen, 5
- AFHCAN software
 - data sharing, 4, 5, 21, 27
 - function, 14, 27
 - interaction with Midmark software, 1, 4
 - limits on backwards compatibility, 27
- automatic diagnosis
 - responsibility to confirm, 2
- Back button, 7
 - on Spirometer screen (image), 5
- battery
 - compartment (image), 2
 - not used on USB models, 37
- Best selection box
 - function and use, 12
 - on data acquisition screen (image), 11
- Bronchodilator pull-down menu
 - function and use, 18
 - on data display screen (image), 17
- Cal button
 - function and use, 13
 - general reference, 29
 - on data acquisition screen (image), 11
- Calibrate button, 6
 - general reference, 19, 29
 - on Spirometer screen (image), 5
- calibration
 - changes in elevation, 21
 - expiration interval, 19
 - function, 19, 29
 - quick verification, 14
 - recommended frequency, 6, 19, 29
 - results out of tolerance, 37
- calibration screen, 7
- calibration syringe
 - adapter, 19
 - function, 4
 - general reference, 2, 35
 - image, 3
 - performing calibration (image), 30
 - procedure, 14, 21
 - type and capacity, 19
- Calibration tab (image), 26
- cleaning, 37
- Configuration tab (image), 24
- Continue Report button, 6
 - general reference, 6, 31
 - on Spirometer screen (image), 6
- Cover Page Settings tab (image), 24
- data acquisition screen, 7
 - image, 11
- data display screen, 7
 - numerical data (image), 17
 - performance curve (image), 17
- Delete Rejected Tests selection box
 - function and use, 13
 - on data acquisition screen (image), 11
- digital spirometer
 - functional description, 1
 - main components, 2
- disposable pneumotach mouthpiece. *See* mouthpiece
- DPM. *See* disposable pneumotach mouthpiece
- Ethnic Adjustments tab (image), 25
- Exit button
 - function and use, 14
 - on calibration screen (image), 21
 - on data display screen (image), 17
- File button
 - function and use, 14
 - on data display screen (image), 17
- flow volume loop curve
 - example (image), 9
- forced vital capacity test. *See* FVC test
- functional description. *See* (device) functional description
- FVC button
 - function and use, 11
 - on data acquisition screen (image), 11
- FVC test
 - description, 1
 - displaying results, 17, 33
 - general reference, 8, 15, 31
 - typical results curve (image), 8
 - viewing multiple results, 9
- handle. *See* spirometer handle
- Help button
 - function and use, 14
 - on calibration screen (image), 21
 - on data acquisition screen (image), 11
 - on data display screen (image), 17
- Incentive/Miscellaneous tab (image), 25
- Interpretation tab (image), 26
- laminar flow element
 - function, 3
 - image, 3

- LED indicator
 - image, 2
 - showing green (significance), 3
- maximum voluntary ventilation test. *See* MVV test
- Measurements tab (image), 26
- Midmark Digital Spirometer. *See* digital spirometer
- Midmark software
 - automatic diagnosis, 2
 - calibration data, 19
 - calibration screen (image), 20, 21
 - data acquisition screen (image), 11
 - data display (image), 17
 - description of on-screen buttons, 13
 - descriptions on-screen buttons, 11
 - elevation and pressure, 20
 - functions, 7
 - general reference, 1, 4, 6, 22, 27
 - interaction with AFHCAN software, 1, 4
 - patient incentive display, 14
 - patient incentive displays (image), 15
 - three major screens (image), 7
- mouthpiece
 - calibration, 19, 21
 - described, 3
 - details (image), 3
 - general reference, 2, 37
 - installed (image), 2
 - laminar flow element. *See* laminar flow element
 - orientation, 4
 - part number, 37
 - precautions, 27
 - proper use, 31
- MVV button
 - on data acquisition screen (image), 11
- MVV test
 - description, 10
 - results curve (image), 10
- New Cal button
 - function and use, 14
 - general reference, 29
 - on calibration screen (image), 21
- Next button
 - function and use, 12
 - on data acquisition screen (image), 11
- nose clip
 - general reference, 2
 - image, 3
 - procedure, 31, 35
- Patient button
 - general reference, 29
- patient feedback message window (image), 16
- patient incentive display image, 15
- Perform Post BD button
 - function and use, 18
 - general reference, 33
 - on data display screen (image), 17
- Perform Pre-BD button
 - function and use, 18
 - general reference, 33
 - on data display screen (image), 17
- Performed By data entry field
 - on calibration screen (image), 21
- Post button
 - function and use, 12
 - on data acquisition screen (image), 11
- power requirements
 - supplied by USB port, 37
- Pre button
 - function and use, 12
 - on data acquisition screen (image), 11
- Pressure data entry field
 - on calibration screen (image), 21
- Prev button
 - function and use, 12
 - on data acquisition screen (image), 11
- Print button
 - function and use, 14
 - general reference, 33
 - on calibration screen (image), 21
 - on data display screen (image), 17
- Reject radio button
 - function and use, 12
 - on data acquisition screen (image), 11
- Reports tab (image), 24
- Review Date data entry field
 - function and use, 19
- Reviewed By data entry field
 - function and use, 19
- Save button, 7
 - general reference, 32, 33
 - on Spirometer screen (image), 5
- Save Review button
 - function and use, 14
 - general reference, 33
- Sensor SN data entry field
 - on calibration screen (image), 20
- Settings button
 - function and use, 13, 18
 - on calibration screen (image), 21
 - on data acquisition screen (image), 11
 - on data display screen (image), 17
- settings dialog box
 - available pages (image), 23
 - description, 23
- spirometer. *See* digital spirometer
- Spirometer button, 5
 - general reference, 5, 29
- spirometer handle
 - cleaning, 37
 - description, 3

- function, 3
- general reference, 2, 19, 31
- image, 2
- Spirometer screen
 - activating Midmark software, 11
 - description of on-screen buttons, 6
 - general reference, 33
 - image, 5
 - initiating calibration, 19, 29
 - report thumbnails, 13
 - returning to, 14, 30
 - supplemental patient information (image), 6
 - viewing thumbnail, 16
- spirometry
 - available tests, 7
 - defined, 1
- Start Cal button
 - function and use, 14
 - general reference, 21, 29
 - on calibration screen (image), 21
- Start New Report button
 - function, 6
 - general reference, 5, 11, 31
 - on Spirometer screen (image), 5
- Start New Test button
 - function and use, 12
 - general reference, 14, 31
 - on data acquisition screen (image), 11
- Stop button
 - general reference, 15
- Syringe S/N data entry field
 - on calibration screen (image), 21
- Syringe Volume data entry field
 - on calibration screen (image), 21
- Test Data button
 - function and use, 13
 - on data acquisition screen (image), 11
- Trending tab (image), 25
- troubleshooting, 37
- VC button
 - on data acquisition screen (image), 11
- VC test
 - description, 9
 - results curve (image), 10
- Verify button
 - function and use, 14
 - on calibration screen (image), 21
- vital capacity test. *See* VC test

