

EBA 20S



Repair instructions

Andreas Hettich GmbH & Co. KG
Föhrenstraße 12, D-78532 Tuttlingen / Germany
Phone (07461) 705-0
Fax (07461) 705-125
info@hettichlab.com, service@hettichlab.com
www.hettichlab.com

© 2007 by Andreas Hettich GmbH & Co. KG

All rights reserved. No part of this publication may be reproduced without the prior written permission of the copyright owner.

Änderungen vorbehalten!

AR2080EN / 03.08

Contents

1	Introduction	5
2	Symbol meanings.....	5
3	Description of the centrifuge.....	6
3.1	Control panel A3	7
3.2	Power board A1	7
3.3	Motor M1	8
3.4	Speed sensor B3.....	8
3.5	Brake resistor R1 and Overtemperature fuse F3.....	8
3.6	Lid lock Y1	8
3.7	Imbalance switch S 2	8
4	Troubleshooting procedures.....	9
5	Error messages.....	10
5.1	MAINS RESET.....	10
5.2	Brief description	10
5.3	Description and elimination of errors.....	11
5.4	Defects without Error indications.....	16
6	Settings and enquiries.....	18
6.1	Setting the machine version.....	18
6.2	Enquiry the machine version.....	19
6.3	Enquiry the programme versions	20
6.4	Brake adjustment	20
6.5	Lid lock control	21
6.6	Imbalance switch-off	21
7	Change mains input fuse.....	22
8	Functional check after a repair.....	22
9	Assembling and disassembling components.....	23
9.1	Motor M1 / speed sensor B3 / Rubber-metal bearing.....	24
9.2	Power board A1	24
9.3	Control panel A3	24
9.4	Brake resistor R1	25
9.5	Overtemperature fuse F3.....	25
9.6	Flat ribbon cable.....	25
9.7	Appliance plug with mains switch.....	25
9.8	Imbalance switch S2	25
9.9	Lid lock Y1	25

10	Circuit diagrams	26
10.1	Circuit diagram and plug assignment.....	27
10.2	Signals and test results.....	28
10.2.1	Speed sensor B3 (Tacho).....	28
10.2.2	Motor	28
10.2.3	Brake resistor R1 and overtemperature fuse F3.....	29
10.2.4	Lid lock Y1	29
10.2.5	Imbalance switch	30
10.2.6	Intermediate voltage on power board A1	30
11	Technical specification	31
12	Skeleton construction of EBA 20S.....	32
13	Index	33

1 Introduction

- Repairs must only be carried out by personnel authorised to do so by the manufacturer.



Interventions and modifications at centrifuges, which have been conducted by persons not authorized by the Andreas Hettich GmbH & Co. KG company, are at their own risk and entail the loss of all guarantee and liability claims. In such an event any guarantee claim or liability claim against the Andreas Hettich GmbH & Co. KG company expires.

- Only original spare parts and original accessories licensed by the Andreas Hettich GmbH & Co. KG company are allowed to be utilised.



If no original spare parts or no original accessories are used, any guarantee claim or liability claim against the Andreas Hettich GmbH & Co. KG company ceases to exist.

- Information about the operation of the centrifuge please see operating instructions.
- We reserve all rights for these technical documents.

2 Symbol meanings



Symbol on the machine:

Attention, general hazard area.

Before using the centrifuge implicitly read the operating instructions and pay attention to the safety relevant references!



Symbol in this document:

Attention, general hazard area.

This symbol refers to safety relevant warnings and indicates possibly dangerous situations.

The non-adherence to these warnings can lead to material damage and injury to personal.



Symbol in this document:

Warning! Danger for human lives by electric shock.



Symbol in this document:

This symbol refers to important circumstances.



Symbol on the machine and in this document:

Symbol for the separate collection of electric and electronic devices according to the guideline 2002/96/EG (WEEE). The device belongs to Group 8 (medical devices).

Applies in the countries of the European Union, as well as in Norway and Switzerland.


3 Description of the centrifuge

The EBA 20S is a microprocessor-controlled centrifuge which is comprised of the following electrical components:

- Control panel
- Power board
- Motor
- Speed sensor (Tacho)
- Brake resistor with overtemperature fuse
- Lid lock system
- Imbalance switch



All electronic components are on mains, due to the DC-coupling

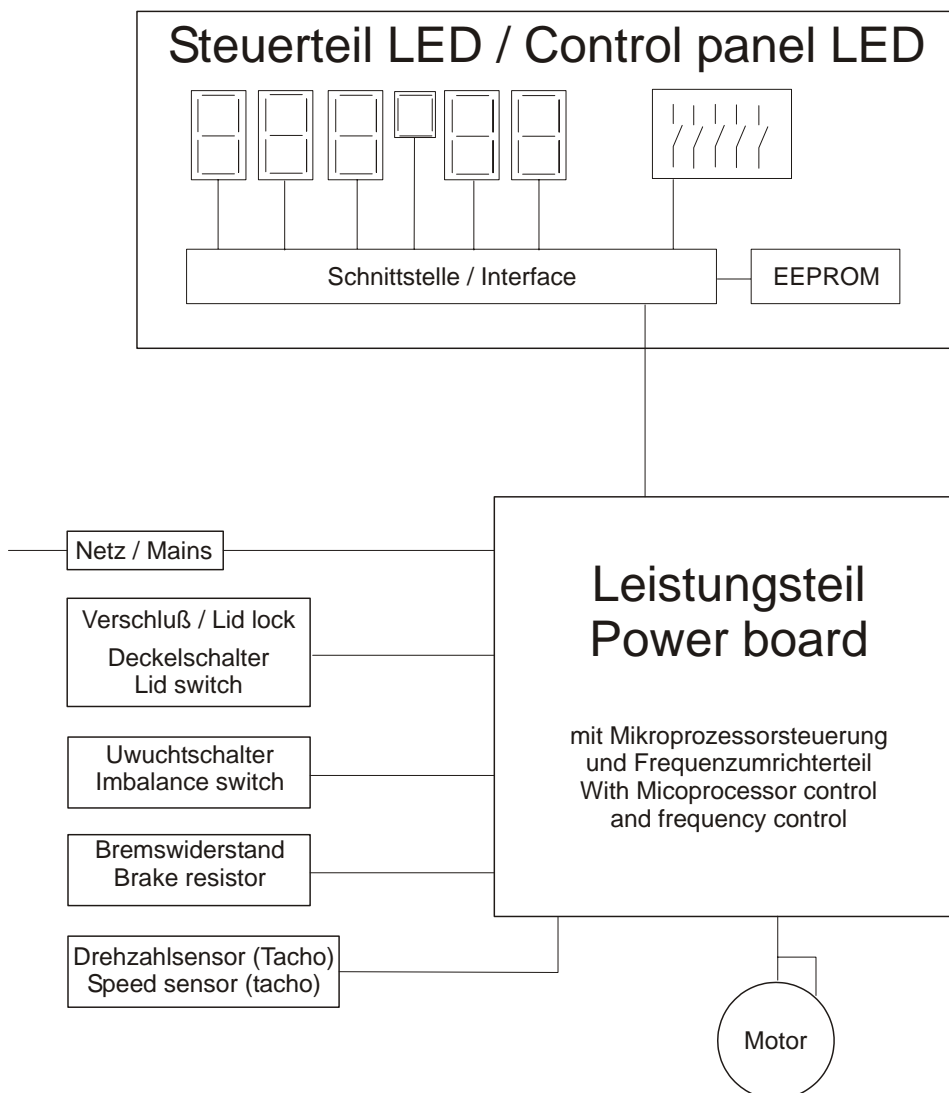

3.1 Control panel A3

The control panel have only restricted control tasks, it disposes of the following characteristics:

- Input panel for operation parameters
- Indication elements
- Transmission of the signals to the power board via the interface.
- Storing of centrifuge version
- Storing the machine version and the brake setting.
By means of the machine version the power board is informed which kind of centrifuge has to be controlled. Then the power board takes the corresponding values from the ROM.
e.g. Max. Speed
Acceleration and deceleration ramps
- Communication with the power board via TTL interface.

The power supply for the control panel is transmitted via the flat ribbon cable:

Pin 1 GND

Pin 4 +5V

3.2 Power board A1

The power board is an combination of:

- Control panel
- Voltage supply
- Frequency converter

The power board carries out the following tasks:

- Power supply 15 V, DC for imbalance switch
- Power supply 5 V, DC for control panel
- Generating the motor power supply.
(three-phase current with variable frequency and voltage)

Functional description: The mains supply is rectified, smoothed and chopped into a pulse width pattern in three bridge elements with a microprocessor.

- Slot for motor with integrated overtemperature switch
- Monitoring the motor current
- Evaluating the overtemperature switch in the motor
- Slot for imbalance switch
- Evaluating the imbalance switch
- Slot for brake resistor
- Triggering the brake resistor
- Slot for speed sensor
- Evaluating the speed sensor pulses (2 pulses per revolution)
- Slot for lid lock
- Triggering the lid lock magnet at standstill of the rotor
- Evaluating the message line lid lock open/closed
- Communication with the control panel via TTL interface
- Error evaluation

3.3 Motor M1

- The motor is a three- phase asynchronous motor with two pairs of poles.
- The motor is protected against overheating by an overtemperature switch.
- The power board evaluates the overtemperature switch.
- The motor is controlled by the power board with a three-phase current with variable frequency and voltage.


3.4 Speed sensor B3

- The speed sensor (speedometer) is located at the bottom of the motor.
- The speed signal (2 pulses per revolution) will be triggered by two magnets fixed at the motor axle.
- The speed of the rotor is monitored and controlled by the power board.

3.5 Brake resistor R1 and Overtemperature fuse F3

- The braking copper which is integrated on the power board transfers the electrical energy produced during braking, from a voltage of 380 V with the 230 V version and 203 V with the 120 V version, to the brake resistor in a controlled manner.
- An overtemperature fuse protects the brake resistor against overheating. When the overtemperature fuse (F3) blows, the power board will be separated from the power supply.

3.6 Lid lock Y1

- Opening of the lid lock is prevented by a latch. The lid lock can only be opened when the relay REL 602 on the power board is energised. This occurs when the rotor is at standstill, mains power is applied and the key  is pressed. The solenoid is energized for three seconds and releases the latch.
- The centrifuge can only be started when the lid is closed. A microswitch on the lid lock detects the position of the lid lock (open/closed) and report it to the power board.

3.7 Imbalance switch S 2


- A switch detects any imbalance.
- Imbalance can only be detected in running mode (starting, centrifuging and braking).
- If any imbalance is detected, the drive is changed over to braking.

4 Troubleshooting procedures

- Fuses in installation in which centrifuge is installed are intact.
- Mains input fuses of centrifuge are intact.
- Supply voltage present at (see circuit diagram):
 - Connecting cable
 - Appliance plug
 - Mains switch
 - Power board A1, plug S402L and S402N
- Look for the displayed error code in the chapter "Error messages".
- Remedy the error according to the instructions.
- Carry out a functional check after every repair and whenever a component is replaced, see chapter "Functional check after a repair".

5 Error messages

The error message will be indicated in the speed display of the front panel
e.g. :

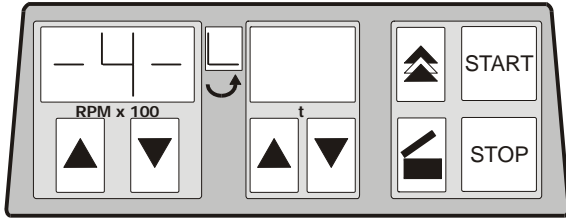

5.1 MAINS RESET

- Switch off the mains switch.
- Wait at least for 10 seconds and then switch on the mains switch again.

5.2 Brief description

Display.	Fault	Brief description	Page
- 1 -	Tacho error	Tacho pulses break down during the run	11
- 2 -	System reset	Mains interrupt	11
- 3 -	Imbalance	Imbalance on the motor axis	11
- 4 -	Communication	Communication error	12
- 5 -	Overload	Fault in the motor or the motor control	12
- 6 -	Overvoltage	Mains is out of the tolerance	13
- 7 -	Overspeed	Overspeed detected	13
- 8 -	Undervoltage	Mains is out of the tolerance	13
- 9 -	Overtemperature	Overtemperature switch in the motor release	14
	Versions Error	No indication in the speed display Wrong machine version setting	16
- c -	Controller-Watchdog	Fault in power board	14
- d -	Lid lock error	Fault in lid lock system	15
- E -	Short circuit	Short circuit in power board	15
- F -	No Rotor code	Wrong machine version	15

5.3 Description and elimination of errors

– 1 – Tacho error

Error	Tacho pulses break down during the run.
Error consequence	Drive switch off and brakes with the adjusted brake level.
Error cause	<ul style="list-style-type: none">• Speed sensor (tacho) defective• Power board (A1) defective• Loosen contact in plug S502
Measurement	<ol style="list-style-type: none">1. A1 / S502 pin 2 (GND) to pin 3 (+U_B)2. Speed sensor plug S502 , pin 4 – pin 2 GND (2 pulses per revolution). See also section 10.2
Error code reset	Wait for a time duration of 120 sec. and after this perform a MAINS RESET.

– 2 – System reset

Error	Mains interrupt during a run
Error consequence	Drive switch off and brakes with the adjusted brake level.
Error cause	<ul style="list-style-type: none">• Power supply has failed• Loosen contact in electrical connections
Error code reset	<ol style="list-style-type: none">1. Wait for rotor standstill2. Open the lid and press key START or perform a MAINS RESET

– 3 – Imbalance

Error	Imbalance on motor axle
Error consequence	Drive switch off and brakes with the adjusted brake level
Error cause	<ul style="list-style-type: none">• Weight difference in rotor components.• Imbalance switch is defective or not connected.• Imbalance switch disadjusted (Adjustment see chapter "Imbalance switch-off").• Loose contact in cable or plug S503• Power board (A1) is defective.
Measurement	Plug of imbalance switch S503 Pin 1 to Pin 4 Switch is a opener. See also chapter 10.2.5.
Error code reset	Open the lid after standstill

– 4 – Communication

Error	Communication error between control board and power board
Error consequence	Drive switch off . No brake effective
Error cause	<ul style="list-style-type: none">• Loose contact in flat ribbon cable• Control board (A3) defective• Power board (A1) defective
Error code reset	MAINS RESET after standstill. If the MAINS RESET performed before standstill, the target speed will be indicated in the display until standstill.

– 5 – Overload

Error	Power board detects overload
Error consequence	Drive switches off. No brake effective
Error cause	<ul style="list-style-type: none">• Power board (A1) defective• Motor defective (to low impedance or motor bearings are defective)
Measurement	Motor coil, see section 10.2.2.
Error code reset	Perform a MAINS RESET after standstill

– 6 – Overvoltage

Error	Overvoltage in intermediate circuit. This error normally only occur when the drive is being braked
Error consequence	Drive switches off. No brake effective
Error cause	<ul style="list-style-type: none">• Power board (A1) defective• Brake resistor defective
Measurement	Measure intermediate voltage, see section 0
Error code reset	Perform a MAINS RESET after standstill.

– 7 – Overspeed

Error	Overspeed. The speed measured by the speed sensor B3 is 250 RPM higher than the maximum speed of the rotor.
Error consequence	Drive switches off. No brake effective
Error cause	<ul style="list-style-type: none">• Speed sensor (B3) defective• Power board (A1) defective
Measurement	Check the speed sensor, see section 10.2.1
Error code reset	Perform a MAINS RESET after standstill.

– 8 – Undervoltage

Error	Mains voltage less than 20% as nominal voltage.
Error consequence	Drive switches off. No brake effective
Error cause	<ul style="list-style-type: none">• Mains voltage too less• Power board (A1) defective
Measurement	Check the mains voltage Measure intermediate voltage, see chapter 0
Error code reset	Perform a MAINS RESET after standstill.

– 9 – Overtemperature

Error	Overtemperature in the motor indicated.
Error consequence	Drive switches off. No brake effective
Error cause	<ul style="list-style-type: none">• Motor defective• Power board (A1) defective• Overtemperature switch in motor defective• Loosen contact in plug S401 Pin 4 or 5
Measurement	Remove plug S405 and measure between Pin 4 and 5 : Switch open : $\approx 0 \Omega$ OK Switch closed: $\infty \Omega$ defective See also section 10.2.
Error code reset	Perform a MAINS RESET after standstill.

– c – Controller-Watchdog

Error	Watchdog in power board Discrepancy in program procedure
Error consequence	Drive switches off. No brake effective
Error cause	<ul style="list-style-type: none">• Power board (A1) defective
Error code reset	Perform a MAINS RESET after standstill.

– d – Lid lock error

Error Lid lock is open during centrifugation.

Error consequence Drive switches off. No brake effective

- Error cause
- Micro switch is defective.
 - Power board (A1) defective.
 - Loosen contact in Plug S404 Pin 5 and 6.
 - Mechanical defect at the lid lock
 - Emergency release active during run.

Measurement Remove plug S404 and measure between Pin 5 and 6 :

Lid closed: $\approx 0 \Omega$

Lid open: $\infty \Omega$

See also section 10.2.

Error code reset Perform a MAINS RESET after standstill.

– E – Short circuit

Error Power board detects short circuit

Current consumption too high

Error consequence Drive switches off. No brake effective

- Error cause
- Power board (A1) defective

Error code reset Perform a MAINS RESET after standstill.

– F – Rotor code

Error Wrong machine version is set.

Error consequence Drive switches off after 15 sec. No brake effective

- Error cause
- Wrong machine version is set.

Error code reset Set the correct machine version, see chapter 6.1


5.4 Defects without Error indications

No speed indication / Machine-Version-Error

Error After switching on the centrifuge the speed indicator extinguish and in the time indicator appears the set machine version.

Error cause • Wrong machine version adjusted

Example for indication:

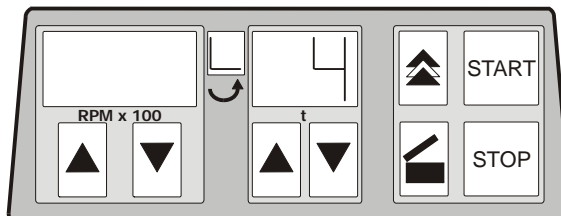

Proceed as the following:

1. Set the machine version 3 (EBA 20S) by using the ▲ and ▼ keys beneath the time indicator.
2. Press the key (STOP) in order to store the set machine version.
3. Perform a MAINS RESET.

The lid can not be opened

Error The lid can not be opened.
The rotation display indicates "Lid opened", but the lid is closed.

Error cause

- Micro switch at lid lock is defective
- Loosen contact on plug S404 / Pin 5 and 6
- Power board A1 defective
- Jumper for lid lock control is not set on power board (A1), see chapter 0

Measurement Remove plug S404 and measure between Pin 5 and 6:

Switch closed: $\approx 0 \Omega$

Switch opened: $\infty \Omega$

See section 10.2.4.

Opening the lid Release the lid by using the emergency release.


The lid can not be opened

Error	The lid can not be opened The rotation display indicates "Lid closed" and the lid is closed.
Error cause	<ul style="list-style-type: none"> • Solenoid on lid lock is defective • Loosen contact on plug S404 / Pin 1 and 3 • Power board A1 defective • Jumper for lid lock control is set on power board (A1), see chapter 0
Measurement	Remove plug S404 and measure between Pin 1 and 3 : Solenoid ok: $\approx 5.7 \text{ k}\Omega$ Solenoid defective: $\infty \Omega$ See also section 10.2.4.
Opening the lid	Release the lid by using the emergency release.

No display

Error	No power supply on control board, no display
Error consequence	No operation possible
Error cause	<ul style="list-style-type: none"> • No mains supply • Power board (A1) defective • Control panel (A3) defective • Flat ribbon cable to control panel (A3) defective • Overtemperature fuse F3 defective
Measurement	<ol style="list-style-type: none"> 1. Mains supply power board A1 S402L to S402N 2. Power board A1 S501 Pin 1 (GND) to S501 Pin 4 (+5V) 3. Control panel A3 S502 Pin 1 (GND) to S502 Pin 4 (+5V) 4. Overtemperature fuse F3. Remove plug S405 and measure directly at both flat connections of the overtemperature fuse: Fuse OK : $\approx 0 \Omega$ Fuse defective: $\infty \Omega$. See also section 10.2.3.

6 Settings and enquiries

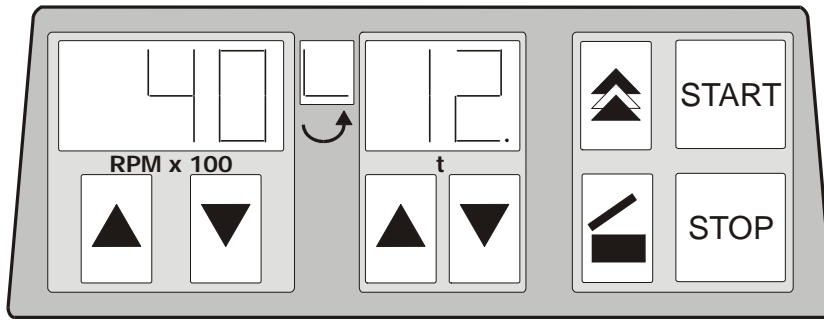

6.1 Setting the machine version

The machine version stored in the control panel must correspond to the centrifuge model. Control panels supplied as spare part are set to machine version 4 (ROTOFIX 32).

After installing the control panel the machine version 3 (EBA 20S) must be set.


If a false machine version is set in the control panel so either after switching on the centrifuge the set machine version will be displayed or after start of the centrifugation run the error message " F –" will be displayed (see chapter "Error messages").

The machine version 3 (EBA 20S) must be set as follows:




All electronic components are on mains, due to the DC-coupling

1. Switch off the mains switch and disconnect the centrifuge from the mains supply.
2. Remove the rotor.
3. Remove the fastening screws of the upper part of the centrifuge housing and place the upper part carefully beside the centrifuge.
4. Short the two throughplatings on the rear side of the control panel with the short circuit plug E1754 :

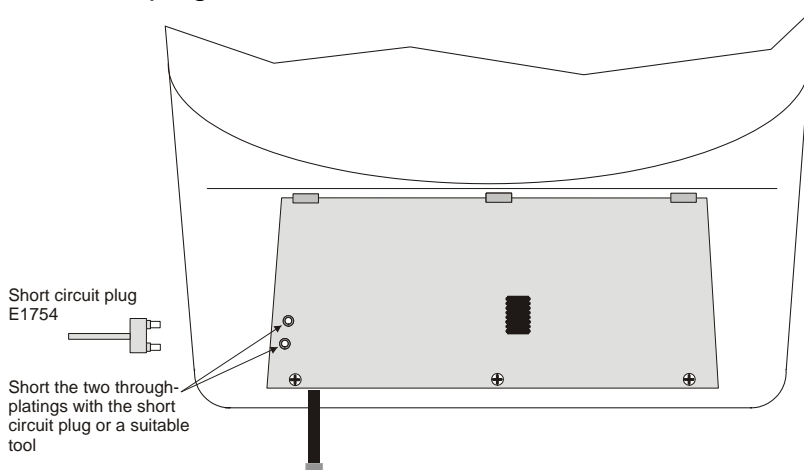


5. Put on carefully the upper centrifuge housing.
6. Connect the centrifuge to the mains and switch the centrifuge on.
7. Switch on the mains switch (position "I").

The set machine version will be displayed, e.g.:

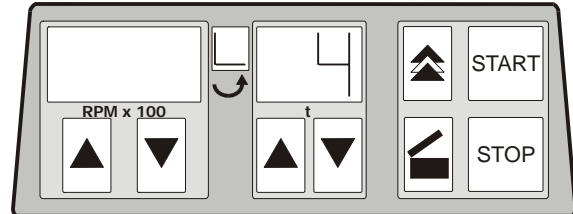


8. Adjust the machine version 3 by using the ▲ and ▼ keys beneath the time indicator.
9. Press the key **STOP** in order to store the set machine version.
10. Switch off the mains switch (position "0").
11. Disconnect the centrifuge from mains and remove the short circuit plug.
12. Put on the upper part of the centrifuge housing and fix it.
13. Switch on the mains switch (position "I").

6.2 Enquiry the machine version

1. Switch off the mains switch.
2. Keep the key ▲ beneath the speed indicator and the key  pressed simultaneously.
3. Switch on the mains switch and release the keys again.

The speed indicator shows the machine version and the time indicator shows the set brake step: e.g.:

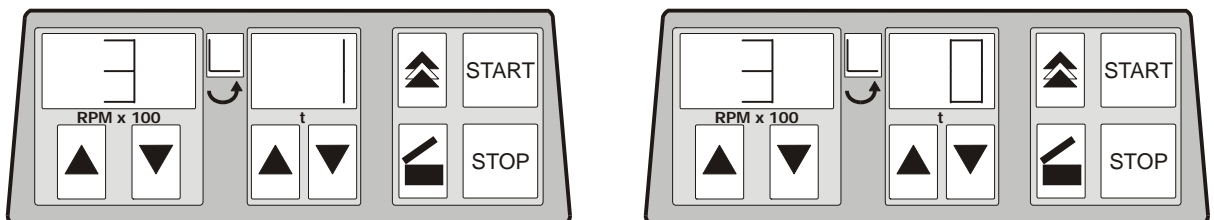


If the machine version and brake step are not displayed, press the ▲ key under the speed indicator until they are displayed.


4. To exit the machine version display press the key **STOP** or perform a MAINS-RESET.

6.3 Enquiry the programme versions

1. Switch off the mains switch.
2. Keep the key ▲ beneath the speed indicator and the key ⏏ pressed simultaneously.
3. Switch on the mains switch and release the keys again.
4. Press the key ▲ beneath the speed indicator so often until the programme version of the control panel (e.g. b1.12) is displayed.

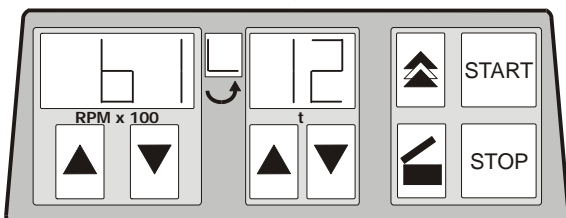



5. Press the key ▲ beneath the speed indicator again. The programme version of the power board (e.g. F1.03) is displayed.

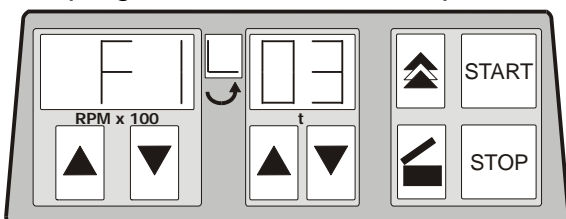

6. To exit the machine version display press the key **STOP** or perform a MAINS-RESET.

6.4 Brake adjustment


1. Switch off the mains switch.
2. Keep the key ▲ beneath the speed indicator and the key ⏏ pressed simultaneously.
3. Switch on the mains switch and release the keys again. The speed indicator shows the machine version and the time indicator shows the set brake step: e.g.:

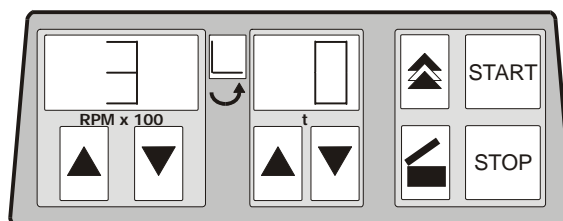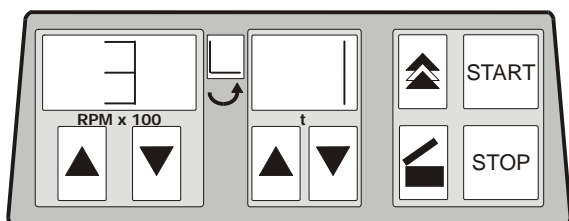


If the machine version and brake step are not displayed, press the ▲ key under the speed indicator until they are displayed.

4. Set the desired brake step with the keys ▲ ▼ beneath the time indicator. Step 1 = short run-down time, Step 0 = long run-down time. For run-down times, see chapter "Anhang/Appendix, Rotor und Zubehör/Rotor and accessories" in the operating instructions.
5. Press the key **STOP** to save the setting.

6.5 Lid lock control

The lid can be opened when the rotor is at standstill and the key  is pressed. The solenoid on the lid lock remains energized for three seconds.

With this kind of lid lock control, the jumper for the lid lock control on the power board (see below) must be plugged.

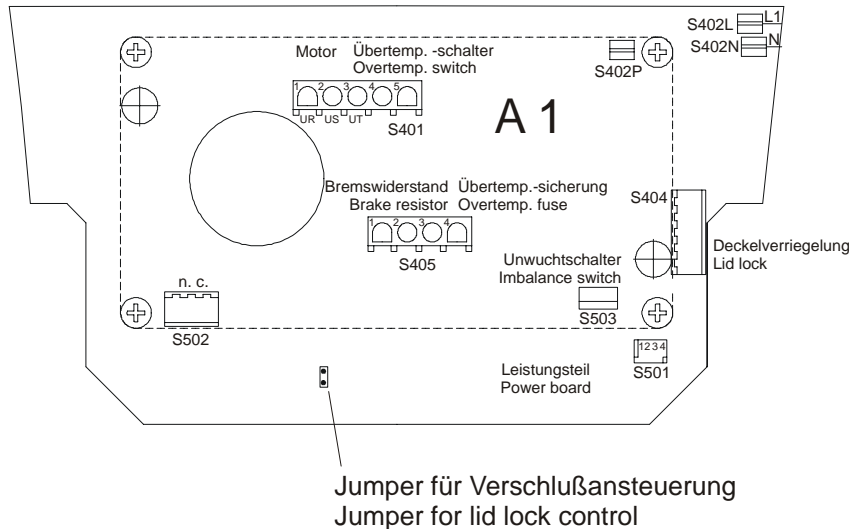

6.6 Imbalance switch-off

The permissible imbalance is specified for rotor E1624 by the indication of the difference in weight of opposite rotor positions.

When having a difference in weight within the range of 7g to 15g during run-up, the drive has to switch off before reaching 1500 RPM.

The imbalance switch-off is adjusted by changing the distance of the imbalance switch.

With a test run with the indicated differences in weight the imbalance switch-off will be checked.


Adjusting the imbalance switch:

- Loosen both screws at the angle bracket of the imbalance switch on the outer part of the housing floor until you can shift it.
- Set the permissible imbalance by shifting the angle bracket.
- Tighten both screws at the angle bracket of the imbalance switch again.
- Check the imbalance switch-off with a test run.

7 Change mains input fuse



Switch off the mains switch and separate the centrifuge from the mains!

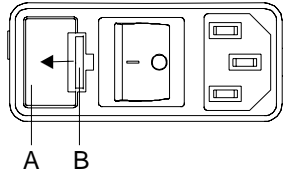

The fuse holder (A) with the mains input fuses is located next to the mains switch.

- Remove the connecting cable from the machine plug socket.
- Press the snap-fit (B) against the fuse holder (A) and remove.
- Exchange defective mains input fuses.



Only use fuses with the rating defined for the type.
See the following table.

- Reinsert the fuse holder until the snap-fit clicks shut.
- Reconnect the centrifuge to the mains supply.

Model	Type	Fuse	Order no.
EBA 20S	2080	T 3,15 AH/250V	E997
EBA 20S	2080-01	T 5 AH/250V	E914

8 Functional check after a repair

After a repair a functional check of the unit must be carried out. For functional check a test run with the loaded rotor must be performed.

During the test run the followings must be checked:

- Function of the keys, the display and the LEDs.
- Run-up and slow-down time, max. speed of the rotor. Values see operating instructions chapter "Anhang/Appendix, Rotoren und Zubehör/Rotors and accessories".
- Sample temperature. Values see operating instructions chapter "Anhang/Appendix, Rotoren und Zubehör/Rotors and accessories".
- Distance of the imbalance switch.
- Current consumption. Values see chapter "Technical specification".

After the test run a safety test must be carried out. Check the following values:

- Insulation resistance > 2 MΩ
 - Protective conductor resistance < 0.2 Ω
 - Leakage current < 3.5 mA *
- * limit according to EN 61010-1

A laboratory centrifuge do not belong to those medical appliances which may be tested according to the regulation IEC 60601-1 or corresponding national medical electronic standards. Laboratory centrifuges are classified as laboratory equipment.

The regulations applying to laboratory equipment are IEC 61010-1 or European standard EN 61010-1.

9 Assembling and disassembling components

Component	Abreversion	Plug	Connected to	Description in section
Motor	M1	S401	A1	9.1
Rubber-metal bearing		--	--	9.1
Power board	A1	--	--	9.2
Control panel	A3	--	--	9.3
Flat ribbon cable	--	--	A1 to A3	9.6
Lid lock	Y1	S404	A1	9.9
Brake resistor	R1	S405	A1	9.4
Overtemperature fuse	F3	S405	A1	9.5
Mains switch	Q1	--	A1	9.7
Appliance plug	--	--	--	9.7
Imbalance switch	S1	S503	A1	9.8



Before assembling or disassembling components, the following working processes must first be carried out to reach the components.

Make a note of the plug numbers !

The components are assembled in reverse order !

- Open the lid
- Switch off the mains switch
- Disconnect the centrifuge from the mains
- Remove the rotor
- Remove the fastening screws of the upper part of the centrifuge housing and place the upper part carefully beside the centrifuge.

9.1 Motor M1 / speed sensor B3 / Rubber-metal bearing



Wait at least 2 minutes after disconnecting the centrifuge from the mains, until the intermediate circuit capacitors of the frequency converter are unloaded.

- Unplug S401 and S502 on the power board.
- Use a socket spanner to loosen and remove the three fastening nuts on lower end plate of motor.
- Lift motor upwards out of centrifuge and unplug the ground wire.
- Unscrew the speed sensor B3 at the bottom of the motor.
- Before the motor is installed, the three rubber-metal bearings must be checked for possible wear or cracks. These items must be replaced if necessary.
- Exchange the motor.
- Care must be taken of the anti-twist device when the rubber metal or motor being installed.

9.2 Power board A1



Wait at least 2 minutes after disconnecting the centrifuge from the mains, until the intermediate circuit capacitors of the frequency converter are unloaded.

- Unplug all plugs and cables from the power board.
- Remove the two fastening screws of the power board
- Exchange the power board.
- Before installation, it must be noted that there is a heat-conducting paste between power board and centrifuge housing floor. Heat conduction from power board to centrifuge housing floor must be ensured.
- The jumper for the lid lock control has to be plugged on the power board (see chapter "Lid lock control").



Heat conduction from power board to centrifuge housing floor must be ensured.

9.3 Control panel A3

- Remove the three fastening screws and lift up the control panel carefully.
- Unplug the flat ribbon cable from the power board.
- Exchange the control panel.

9.4 Brake resistor R1



Wait at least 2 minutes after disconnecting the centrifuge from the mains, until the intermediate circuit capacitors of the frequency converter are unloaded.

- Remove the two fastening screws.
- Unplug the two cable from the brake resistor.
- Exchange the brake resistor.

9.5 Overtemperature fuse F3

- Remove the two cables from the overtemperature fuse.
- Unscrew the fastening screws.
- Exchange the overtemperature fuse.

9.6 Flat ribbon cable

- Unplug the flat ribbon cable from the power board.
- Disassemble the control board (see section 9.3).
- Unplug the flat ribbon cable from the control board.
- Assemble in reverse order.

9.7 Appliance plug with mains switch

- Press the upper and the lower fastening clip and remove the appliance plug out of the housing
- Remove the plugs and the earth connection from the appliance plug.
- Remove the fuse holder, see chapter 7.
- Exchange the appliance plug.

9.8 Imbalance switch S2

- Unplug S503 from the power board.
- Remove the to fastening screws.
- Remove the cables at the imbalance switch.
- Exchange the complete imbalance switch.
- After mounting the imbalance switch, adjust it as described in chapter "Imbalance switch-off".

9.9 Lid lock Y1




- Remove the two fastening screws on the top of the upper housing.
- Unplug the cable from the switch and the solenoid.
- Remove the earth connection from the lid lock.
- Exchange the lid lock.

10 Circuit diagrams

Used cable colours and their abbreviations:

Abbreviation	Colour
BK	black
BN	brown
BU	blue
GN/YE	green-yellow
GY	grey
RD	red
WH	white
YE	yellow

10.1 Circuit diagram and plug assignment

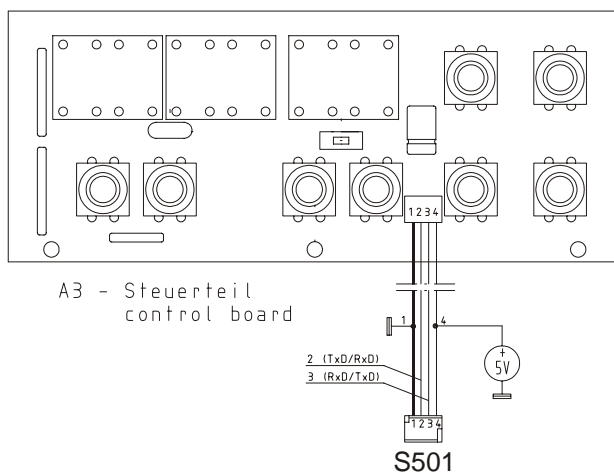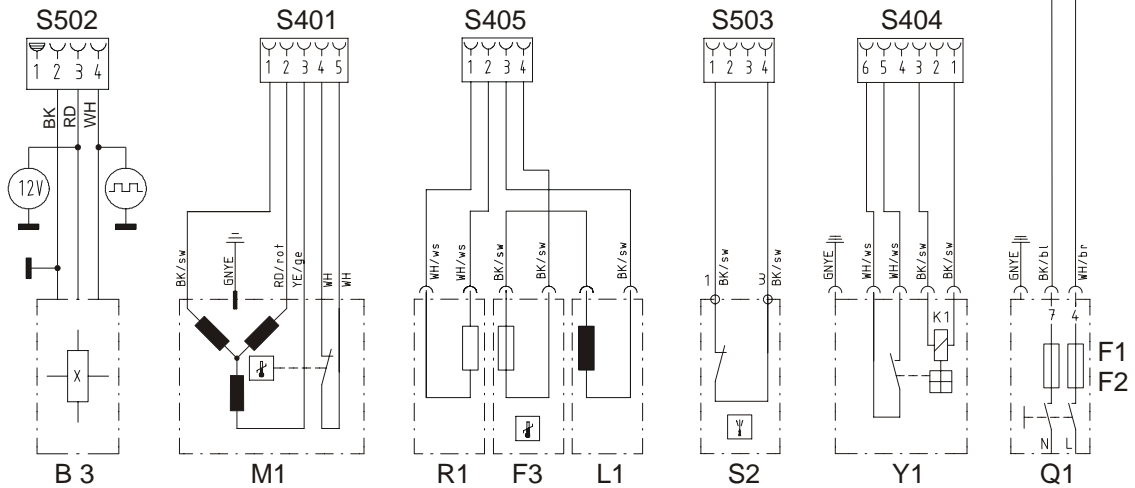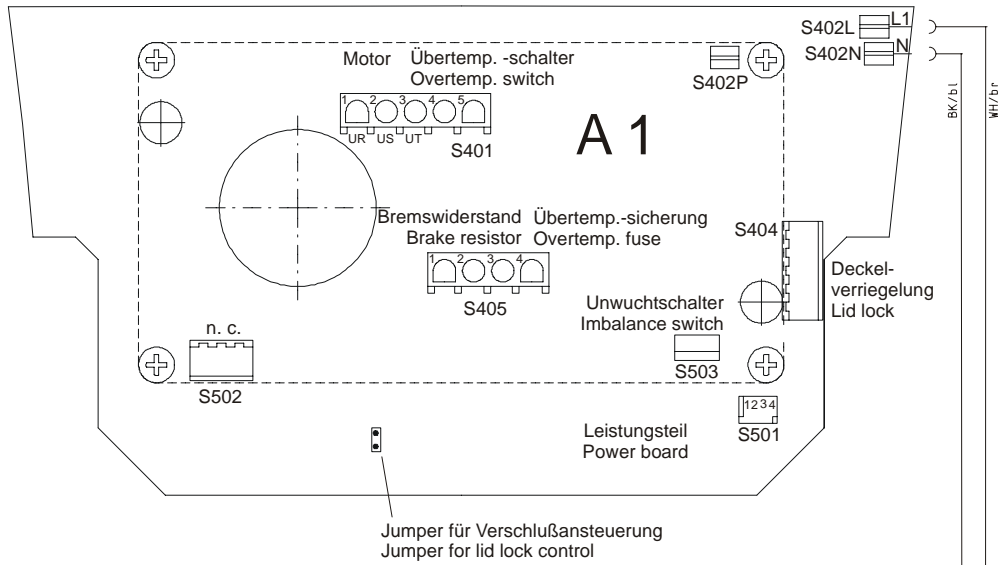

Typ / type	F1, F2
2080	T 3.15 AH/250V
2080-01	T 5 AH/250V

10.2 Signals and test results

10.2.1 Speed sensor B3 (Tacho)

Power supply:

Plug S502, measure between Pin 2 (GND) and Pin 3:

+8...+15V

0V


Tacho signal:

Plug S502 measure between Pin 2 (GND) and Pin 4:

+8...+15V

0V

2 pulses per revolution

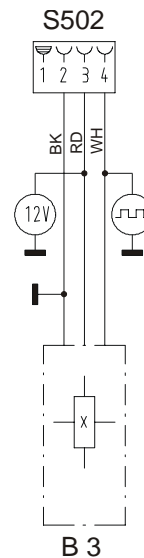

10.2.2 Motor

Motor coil

Remove plug S401 and measure between two lines:

230 V version: $\approx 8.4 \Omega$


120 V version: $\approx 2.1 \Omega$

Overtemperature switch

Remove plug S401 and measure between pin 4 and 5 :

Switch closed : $\approx 0 \Omega$

Switch released: $\infty \Omega$

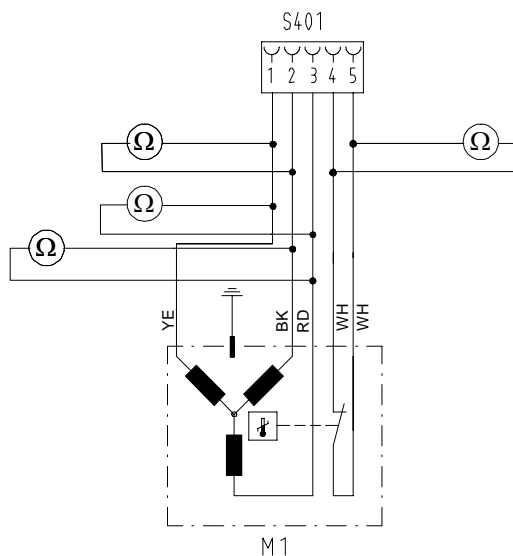

10.2.3 Brake resistor R1 and overtemperature fuse F3

Brake resistor R1

Remove plug S405 and measure between pin 1 and 2 :

230 V Version: $\approx 330 \Omega$


120 V Version: $\approx 82 \Omega$

Overtemperature fuse F3

Remove plug S405 and measure between pin 3 and 4 :

Fuse OK : $\approx 0 \Omega$

Fuse defective: $\infty \Omega$

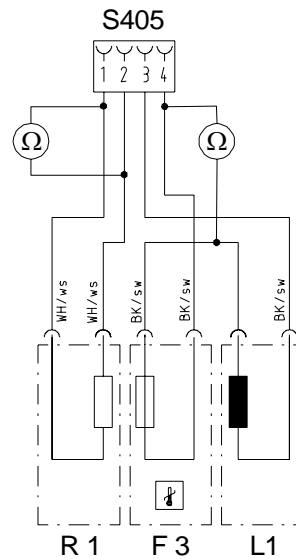

10.2.4 Lid lock Y1

Lid switch:

Remove plug S404 and measure between pin 5 and 6 :


Switch pressed : $\approx 0 \Omega$

Switch inoperative: $\infty \Omega$

Lid lock solenoid K1

Remove plug S404 and measure between pin 1 and 3 :

230 V Version: $\approx 5.7 \text{ K}\Omega$

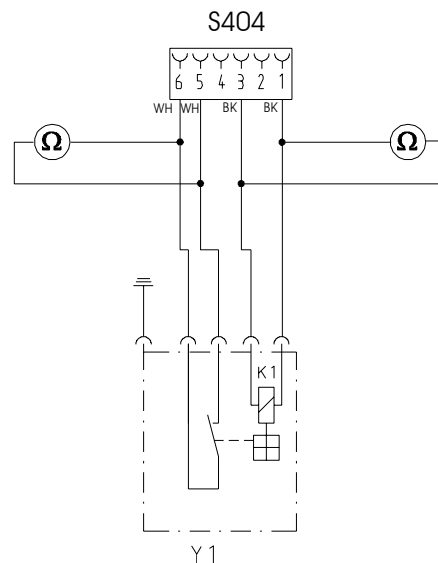


10.2.5 Imbalance switch

Remove plug S503 and measure between pin 1 and pin 4:

Switch pressed : $\infty \Omega$

Switch inoperative: $\approx 0 \Omega$

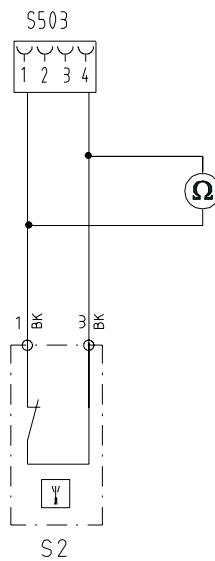

10.2.6 Intermediate voltage on power board A1


The intermediate voltage has to be measured during the run (run-up, run and run-down).

230 V version:

380 V DC

120 V version:


203 V DC

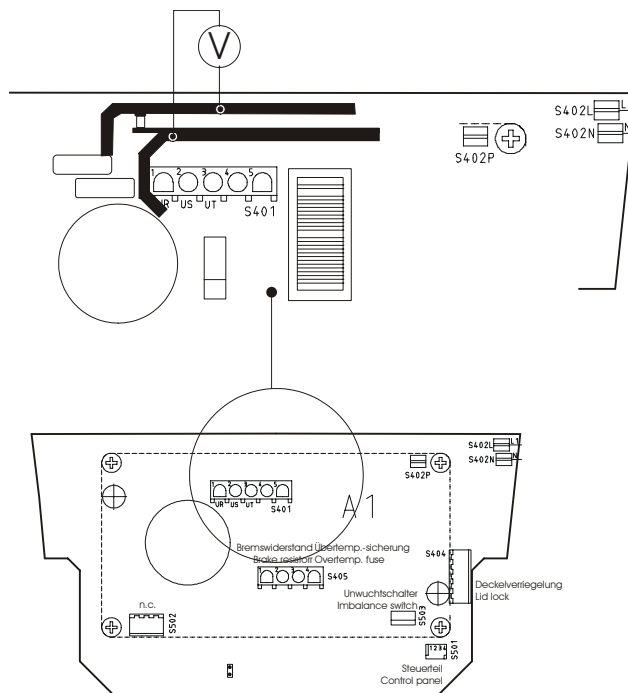

11 Technical specification

Manufacturer	Andreas Hettich GmbH & Co. KG D-78532 Tuttlingen	
Model	EBA 20S	
Type	2080	2080-01
Mains voltage ($\pm 10\%$)	200 - 240 V 1~	100 - 127 V 1~
Mains frequency	50 - 60 Hz	50 - 60 Hz
Connected load	130 VA	180 VA
Current consumption	0.7 A	1.8 A
Max. capacity	8 x 15 ml	
Allowed density	1.2 kg/dm ³	
Speed (RPM)	8000	
Force (RCF)	6153	
Kinetic energy	1450 Nm	
Obligatory inspection (BGR 261)	no	
Ambient conditions (EN 61010-1)	<p style="text-align: center;">Indoors only Up to 2000 m above sea level 2°C to 40°C Maximum relative humidity 80% for temperatures up to 31°C, linearly decreasing to 50% relative humidity at 40°C.</p>	
– Set-up site		
– Altitude		
– Ambient temperature		
– Humidity		
– Excess-voltage category (IEC 60364-4-443)	II	
– Pollution degree	2	
Device protection class	I	
Not suitable for use in explosion-endangered areas.		
EMC		
– Emitted interference (suppression of radio interference)	EN 55011, Group 1, Class B EN 61000-3-2 EN 61000-3-3	FCC Class B
– Interference immunity	EN 61000-6-1	----
Noise level (dependent on rotor)	≤ 55 dB(A)	
Dimensions		
– Width	231 mm	
– Depth	292 mm	
– Height	216 mm	
Weight	approx. 6.3 kg	

12 Skeleton construction of EBA 20S

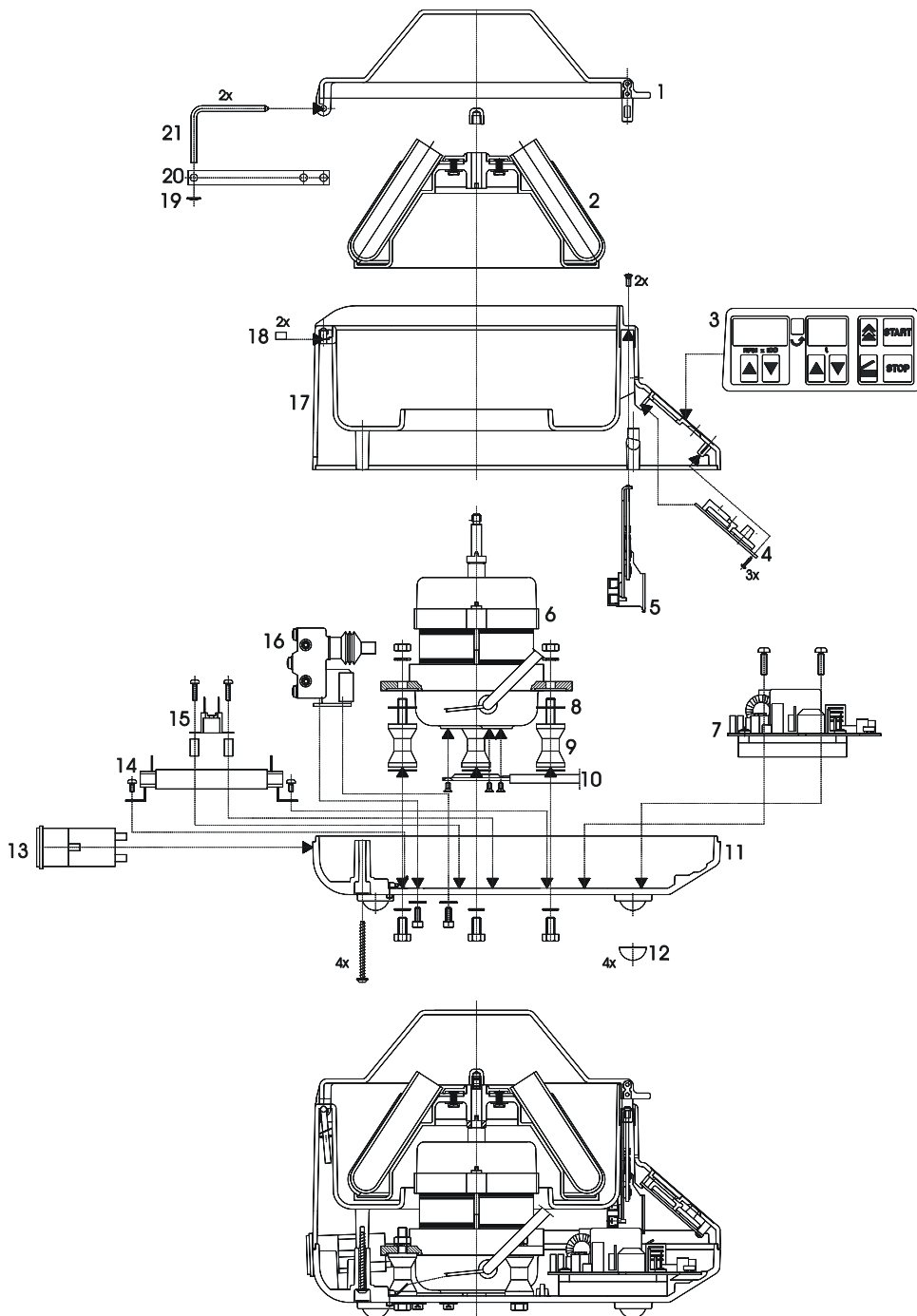

- 1 Lid complete
- 2 Rotor
- 3 Covering foil for control panel
- 4 Control panel
- 5 Lid lock
- 6 Motor
- 7 Power board
- 8 Anti-twist device
- 9 Rubber-metal bearing
- 10 Speed sensor
- 11 Housing lower part
- 12 Rubber foot
- 13 Appliance plug
- 14 Brake resistor
- 15 Overtemperature fuse
- 16 Imbalance switch
- 17 Housing upper part
- 18 Friction rubber
- 19 Lock washer
- 20 Flat spring
- 21 Hinge bolt

13 Index

- 1 – Tacho error	11
- 2 – System reset	11
- 3 – Imbalance	12
- 4 – Communication	12
- 5 – Overload	12
- 6 – Overvoltage	13
- 7 – Overspeed	13
- 8 – Undervoltage	13
- 9 – Overtemperature	14
- c – Controller-Watchdog	14
- d – Lid lock error	15
- E – Short circuit	15
- F – Rotor code	15
Appliance plug	25
Assembling	23
Brake adjustment	20
Brake resistor	8, 25, 29
Cable colours	26
Circuit diagram	27
Control panel	7, 24
Controller-Watchdog – c –	14
Description and elimination of errors	11
Disassembling	23
Enquiries	18
Error messages	10
Flat ribbon cable	25
Functional check	22
Imbalance	12
Imbalance switch	8, 25, 30
Imbalance switch-off	21
Intermediate voltage	30
Lid lock	8
Lid lock control	21
Lid lock solenoid	29
Lid lock Y1	25
Lid switch	29

Machine version	18
Machine-Version-Error	16
Mains input fuses	22
Mains switch	25
Motor	8, 24, 28
Motor coil	28
No display	17
No speed indication	16
Overtemperature fuse	8, 25
Overtemperature switch	28
Plug assignment.....	27
Power board.....	7, 24, 30
Settings	18
Signals	28
Speed sensor	8, 11, 24
Test results	28
The lid can not be opened.....	17