


GRASEBY

omnifuse

Syringe Pump

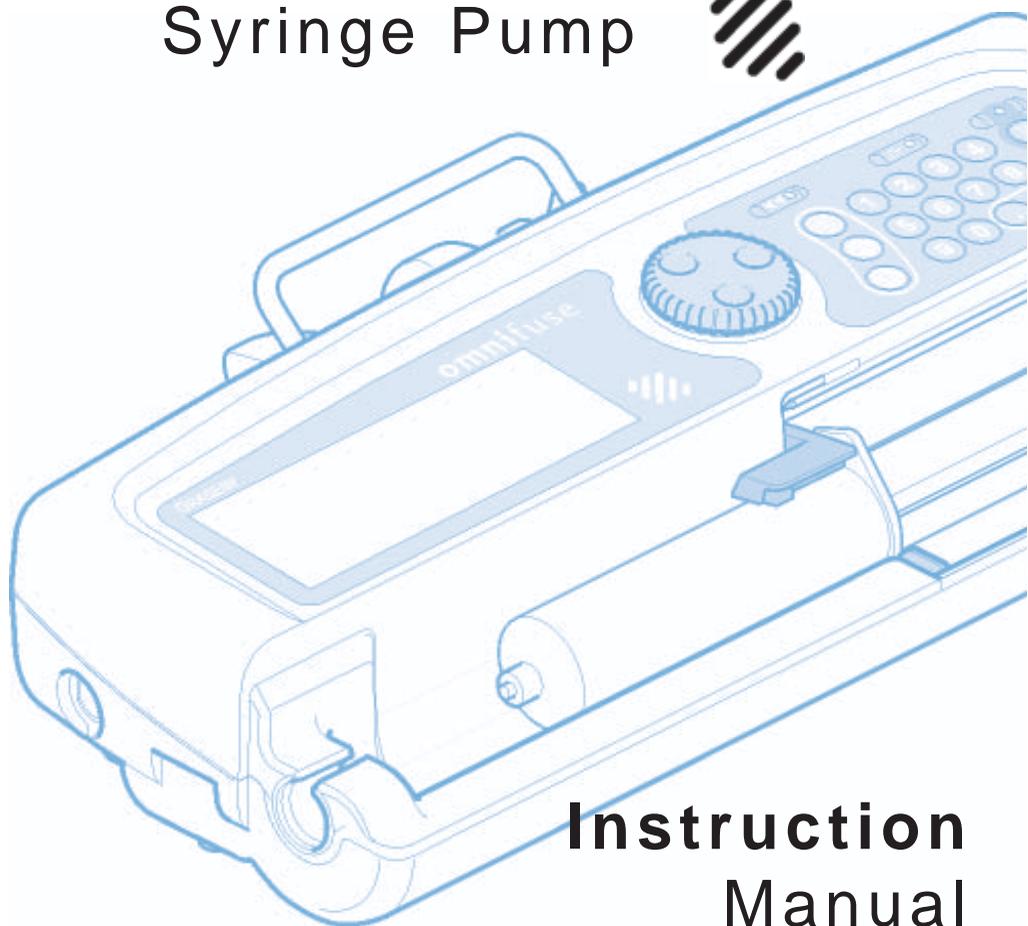

**Instruction
Manual**

Published by Graseby Medical Limited.

All possible care has been taken in the preparation of this publication, but Graseby Medical Limited accepts no liability for any inaccuracies that may be found.

Graseby Medical reserves the right to make changes without notice both to this publication and to the product which it describes.

Copyright © 2004 Graseby Medical Limited

No part of this publication may be reproduced, transmitted, transcribed, or stored in a retrieval system or translated into any human or computer language in any form or by any means without the prior permission of Graseby Medical Limited.

GRASEBY MEDICAL LIMITED

Colonial Way,
Watford,
Hertfordshire,
United Kingdom,
WD24 4LG

Telephone: +44 (0)1923 246434

Facsimile: +44 (0)1923 231595

<http://www.smiths-medical.com>

Registered in England. Company number 995550

Trademarks and acknowledgements:

“Graseby”; “Smiths”; “omnifuse”; “omnistack”;  and “Flo-Safer” are all trademarks of the Smiths Group plc.

All other trademarks are acknowledged as the property of their respective owners.

Using this manual

Read the first three chapters of this manual for an overview of the Omnifuse pump and to find out how to program the pump for a basic infusion.

Warnings and Cautions

This gives a list of Warnings and Cautions which you must be aware of before using the pump. These are repeated on the relevant pages in the manual.

Table of Contents

Use the table of contents to see the structure of the manual and the headings in the order in which they appear.

Chapter 1, Introduction

This introduces the Omnifuse pump. It provides a high-level description of the features of the pump.

Chapter 2, Basics

This chapter covers the external features of the pump, and contains detailed explanations of the techniques for using the pump, for example, switching on and loading or changing a syringe. It also covers the care and maintenance of the pump.

Chapter 3, Programming the pump

This chapter describes how to program the Omnifuse pump for a Continuous infusion. It explains how to use the pump features designed for the clinical environment, for example, setting the occlusion alarm, viewing and resetting totals, purging the line, using Sleep mode, and giving a bolus.

Read the remaining chapters for details on programming and running other types of infusion and to find out how to use the pump variants.

Chapter 4, Specialised infusions

This chapter explains how to program using Mass Units, and describes how to use the Preset Time, Preset Volume, Intermittent and Circadian Rhythm infusion modes on the Omnifuse.

Chapter 5, Protocols

This chapter explains how to program the Omnifuse by selecting from a menu of Drug Protocols.

Chapter 6, In-line occlusion sensing

This chapter explains how to use an Omnifuse pump variant with an In-line occlusion sensor.

Chapter 7, Lockable cover

This chapter explains how to use the Omnifuse pump variant with a Lockable cover.

Chapter 8, Troubleshooting

This chapter explains the warnings and alarms that may be displayed by the Omnifuse.

Specifications and Standards

This lists the specifications for the Omnifuse pump and the standards with which it complies.

Index

At the end of the manual is an index which provides a list of key words, and cross-references these to the relevant pages in the manual.

Warnings/Cautions

Warnings

Warnings tell you about dangerous conditions that could lead to death or serious injury to the user or patient, that can occur if you do not obey all of the instructions in this manual.

1. **WARNING:** To avoid over- or under- infusion, always verify that the brand and size of the loaded syringe are the same as the brand and size displayed on the screen before starting an infusion. Failure to do so may result in an inaccurate delivery of medication, resulting in patient injury or death.
2. **WARNING:** To avoid incorrect or inappropriate configuration of the pump, the Configuration must only be changed by qualified persons or authorised personnel. Incorrect pump configuration could lead to inappropriate infusion resulting in patient injury or death.
3. **WARNING:** This equipment is not suitable for use in the presence of flammable anaesthetic mixtures with air, oxygen, or nitrous oxide. The use of the device in presence of such mixtures may lead to explosion or fire.
4. **WARNING:** To avoid possible malfunction of the pump, do **not** expose the pump to X- rays, gamma rays or ionizing radiation, or to the RF interference or strong electric/magnetic fields emitted (for example) by diathermy equipment or mobile telephones. If the pump is used in the presence of, or in combination with Magnetic Resonance Imaging (MRI) machines it must be protected from the magnetic field emitted by such equipment. Malfunction of the pump can cause incorrect infusion or loss of infusion resulting in patient injury or death.
5. **WARNING:** Operation of the pump outside the temperature limits defined in the specification may result in erroneous operation. Ensure that the temperature is within the specified limits. Failure to do so may result in patient injury or death.
6. **WARNING:** In order to ensure that the intended infusion is performed, data must be entered correctly. Before confirming any displayed data the user should ensure that it is correct. Failure to do so may result in compromised function of the product, patient injury or death.
7. **WARNING:** Failure to follow the Maintenance Procedures described in Chapter 2 of the *Omnifuse Service Manual* may result in compromised function of the product and lead to patient injury or death.
8. **WARNING:** It is essential that clinical staff remain within visual and audible range of the pump so that they can respond promptly to critical alarms. Failure to respond promptly to an alarm may result in patient injury or death.
9. **WARNING:** The user should ensure that the performance offered by the pump is fit for the intended purpose. Failure to do so may result in compromised function of the product, patient injury or death.
10. **WARNING:** When the pump is carrying out an infusion, to ensure that electrical safety is maintained only items of equipment that conform to EN60950 are to be connected to the RS232 connector situated at the base of the pump, otherwise patient safety may be compromised.
11. **WARNING:** While Graseby Medical Limited have taken all reasonable steps to ensure that the pump operates correctly while under remote control, it is the responsibility of the person who designs and implements the controlling device to ensure that the resulting system (pump and controlling device) is fit for its intended purpose. Failure to do so may result in compromised function of the product, patient injury or user injury.
12. **WARNING:** Do not use a faulty pump. If the pump detects a fault, a continuous alarm will sound and the screen will display a System Fault message. If this happens, switch the pump off, disconnect it from the mains and take it to a suitably qualified engineer. Incorrect performance of the pump can cause complications resulting in patient injury or death.
13. **WARNING:** Failure to use the mains lead supplied with the pump will compromise the pump's ability to resist fluid ingress, resulting in possible user or patient injury or death.
14. **WARNING:** Correct management of battery charging is essential to ensure that the pump can operate on batteries for the time specified. Failure to do so may lead to impaired functioning of the pump, resulting in patient injury or death.
15. **WARNING:** The occlusion alarm level must be checked before starting an infusion to ensure that it is appropriate for that infusion. Failure to do so may result in an unacceptably slow time to occlusion alarm, resulting in patient injury or death.

16. **WARNING:** If an occlusion alarm occurs, immediately clamp the line to eliminate the possibility of a bolus being delivered to the patient. Then inspect the fluid pathway for kinks, clogged catheter, etc. in order to remove the occlusion prior to restarting the infusion. An unintentional bolus of medication can result in patient injury or death.
17. **WARNING:** Use only the syringes and administration sets listed in *Specifications and Standards* at the end of this manual. Failure to do so may result in an inaccurate delivery. Graseby Medical does not guarantee performance of the pump if syringes other than those listed are used. Incorrect function or performance of the pump can cause complications resulting in patient injury or death.
18. **WARNING:** When using the In-line Occlusion Sensing option, use only a Graseby Medical or Graseby Medical-approved extension set (part number 0130-0041). Graseby Medical does not guarantee performance of the pump if unapproved extension sets are used. Failure to observe this warning may lead to compromised performance of the pump, resulting in patient injury or death.
19. **WARNING:** The volume of fluid contained in the connecting tubing is a residual amount and will **not** be infused. Allowance must be made for this extra volume of fluid when initially filling the syringe and purging the system. Under-delivery of medication can cause complications resulting in patient injury or death.
20. **WARNING:** To avoid patient embolism, ensure that the patient tubing is purged of all air bubbles before administering any medication. The pump provides a purge facility to assist with this process. The presence of air within the medication can result in complications leading to patient injury or death.
21. **WARNING:** For safe operation of the pump, the syringe must be correctly loaded. Ensure that the syringe plunger is properly aligned before closing the barrel clamp. Failure to do so may result in inaccurate delivery, resulting in patient injury or death.
22. **WARNING:** Ensure that your fingers are not in the path of the pusher during syringe loading or unloading. Failure to do so may result in user injury.
23. **WARNING:** To avoid syphoning of the syringe contents (free flow), ensure that the syringe is correctly loaded into the pump, that the syringe plunger is properly engaged by the pump's pusher block and that the pump is placed not more than 80cm above the infusion site. Syphoning can result in over-infusion leading to patient injury or death.
24. **WARNING:** To avoid over-infusion, do not purge the infusion line when the administration set is connected to the patient. Over-infusion of medication can result in patient injury or death.
25. **WARNING:** To avoid the pump becoming detached from an IV pole always make sure that the pump is securely fixed to the pole. Always check the security and stability of the assembly with the pump mounted. If no IV pole is used make sure that the pump is completely stable on a horizontal surface. Failure to observe this warning may cause damage to the Omnifuse pump and harm the user or the patient. As a result, the user or patient may suffer direct injury, or the Omnifuse pump may fail to operate correctly, leading to patient injury or death.
26. **WARNING:** Following a significant liquid spill onto the pump, it should be wiped dry and inspected by service personnel before being returned to service. Failure to do so may result in compromised functioning of the pump, leading to patient or user injury or death.
27. **WARNING:** To avoid syphoning of the syringe contents (free flow), ensure that the patient line is clamped before loading or unloading the syringe. Syphoning can result in over-infusion leading to patient injury or death.
28. **WARNING:** When the Omnifuse pump is fitted with a Lockable cover, the cover and associated security software should always be used in accordance with local protocol. The Lockable cover protects the syringe and its contents from tampering; the security software ensures that the infusion parameters are not interfered with. Failure to use the security software together with the Lockable cover could result in an inaccurate delivery leading to patient injury or death.

Cautions

Cautions tell you about dangerous conditions that may occur and cause damage to the pump if you do not obey all of the instructions in this manual.

1. **CAUTION:** Refer all service, repair and calibrations only to qualified technical personnel. Unauthorised modifications to the pump must **not** be carried out.
2. **CAUTION:** When turning the pump on, if screens similar to those illustrated are not displayed do not use the pump, and send the pump to authorised service personnel.
3. **CAUTION:** Do not attempt to move the pump's pusher by hand. Always use the syringe Load key  to move the pusher. Failure to observe this caution may cause mechanical damage to the pump.
4. **CAUTION:** Never carry the pump except by the handle. Failure to do so may result in damage to the case, or you may drop the pump and cause it internal damage.
5. **CAUTION:** Do not use cleaning and disinfecting agents other than the approved ones specified here.
6. **CAUTION:** The pump must **not** be immersed in any liquids or exposed to strong organic solvents. Wipe off spills immediately, and do not allow fluid or residues to remain on the pump. Additionally, the pump is not designed to be autoclaved, steam-sterilised, EtO-sterilised or subjected to temperatures in excess of 45° C (113° F). Failure to observe this caution may cause serious damage to the pump.
7. **CAUTION:** Failure to use the mains lead retainer means that the pump may be accidentally disconnected from the AC mains supply.
8. **CAUTION:** Users should bear in mind that the syringe-ear clamp is for location only and may not be powerful enough to hold the syringe in place against the powerful negative back-pressures that may be encountered in certain clinical applications.

Contents

Chapter 1 - Introduction

- Introduction to the Omnifuse pump 1 - 1**
 - Omnifuse pump features 1 - 1
 - Omnifuse pump variants 1 - 4
- PC software 1 - 5**
- Omnistack pump stacking system 1 - 6**

Chapter 2 - Basics

- Parts of the Omnifuse pump 2 - 1**
 - Rear of the pump 2 - 1
 - Front of the pump 2 - 2
 - The Command wheel, keypad and LEDs 2 - 3
- Switching the pump on and off 2 - 4**
 - Switching on and power-up tests 2 - 4
 - Switching off 2 - 5
 - Using the pump on batteries 2 - 6
 - Recharging the batteries 2 - 7
- Omnifuse screen, icons and indicators 2 - 8**
 - Icons on the pump screen 2 - 8
 - Screen icons 2 - 9
 - Pump indicators 2 - 10
 - Sounds on the Omnifuse 2 - 11
 - Setting the alarm volume 2 - 11
- Loading and unloading a syringe 2 - 12**
 - When to load the syringe 2 - 12
 - Syringe brand and size 2 - 13
 - Replacing a syringe during an infusion 2 - 13
 - Warnings and Cautions concerning syringes 2 - 14
 - How to load a syringe 2 - 15
 - How to unload a syringe 2 - 16
 - How to change the brand of syringe 2 - 17
- Care and maintenance 2 - 18**
 - Carrying the pump 2 - 18
 - Positioning the pump for use 2 - 18
 - Using the pole clamp 2 - 19
 - Cleaning and care of the pump 2 - 20

Contents

Chapter 3 - Programming the pump

Continuous infusion in ml per hour	3 - 1
Changing the rate while running	3 - 2
Ending a Continuous infusion	3 - 2
Occlusion pressure and alarms	3 - 3
How to set or change the occlusion alarm level	3 - 4
Monitoring the dry-side pressure	3 - 5
The Totaliser	3 - 6
About the cumulative total	3 - 6
How to view and reset the cumulative total	3 - 7
Hourly total and Shift total	3 - 8
Purging the line	3 - 10
Sleep mode	3 - 11
How to use SLEEP mode	3 - 11
While the pump is asleep	3 - 11
Waking the pump up	3 - 12
Omnifuse bolus feature	3 - 13
Types of bolus delivery methods	3 - 13
Default rate for bolus	3 - 13
Bolus beep	3 - 13
Quitting the Bolus screen	3 - 14
Bolus and TOTAL	3 - 14
How to give a hands-on bolus	3 - 15
How to give a pre-set bolus	3 - 16
Graphics	3 - 17
Accessing the Graphics screen	3 - 17
Viewing infusion data as a graph	3 - 18
Graph icons	3 - 18
How to view an Infusion Profile graph and text	3 - 19

Chapter 4 - Specialised infusions

Mass units and special infusion modes	4 - 1
Setting up the mass units and extra modes	4 - 1
Programming using mass units	4 - 2
Setting up a Continuous infusion in mg/h	4 - 3
Setting up a Continuous infusion in $\mu\text{g}/\text{kg}/\text{min}$	4 - 4
Setting up a Continuous infusion with user mix drug concentration	4 - 5

How to change the infusion mode 4 - 7

The Review screen 4 - 8

Preset Volume infusion 4 - 9

 Changing the rate while running - Preset Volume 4 - 10

 Checking the totals - Preset Volume 4 - 11

 Giving a bolus - Preset Volume 4 - 11

 Ending a Preset Volume infusion 4 - 11

Preset Time infusion 4 - 12

 Changing the rate while running - Preset Time 4 - 13

 Giving a bolus - Preset Time 4 - 14

 Checking the totals - Preset Time 4 - 14

 Ending a Preset Time infusion 4 - 14

Intermittent infusion 4 - 15

 Programming an Intermittent infusion 4 - 16

 Starting an Intermittent infusion 4 - 17

 Background rates and KVO rate 4 - 18

Circadian Rhythm infusion 4 - 19

 Setting up and running a Circadian Rhythm infusion 4 - 20

Chapter 5 - Drug Protocols

About Omnifuse Drug Protocols 5 - 1

 Programming using a protocol 5 - 1

Selecting a protocol 5 - 2

 Programming the pump without a protocol 5 - 3

..... 5 - 4

Chapter 6 - In-line occlusion sensing

In-line Occlusion Sensing 6 - 1

 In-line pressure 6 - 1

Using the In-line pressure sensor 6 - 2

Monitoring In-line pressure 6 - 3

 How to view a Pressure Trend graph 6 - 4

 In-line sensor warning messages 6 - 6

Contents

Chapter 7 - Lockable cover

Introduction	7 - 1
Features of the Lockable cover	7 - 1
About the security software	7 - 2
Infusions with the cover and security software	7 - 3
Removing and replacing the cover	7 - 5

Chapter 8 - Troubleshooting

Introduction	8 - 1
Format of Omnifuse Alarms and Warnings	8 - 1
How to handle Alarms	8 - 1
How to handle Warnings	8 - 2
How to handle System faults	8 - 3
Alarms	8 - 4
Warnings	8 - 6
One-off warnings	8 - 6
Repeated warnings	8 - 11


Specification

Specification	S - 1
Standards	S - 7
Startup curves	S - 8
Trumpet curves	S - 9

Index

Chapter 1

Introduction

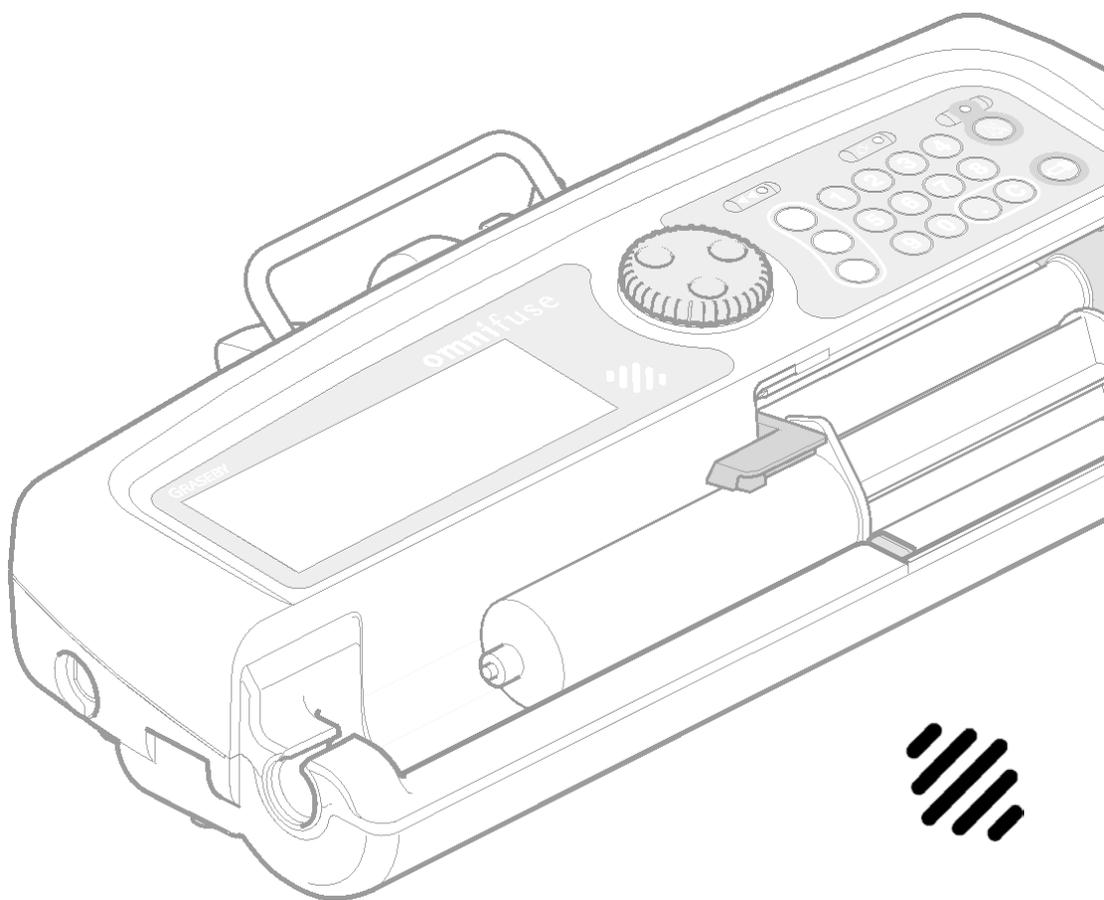

Introduction to the Omnifuse pump

The Omnifuse is part of the wide range of Graseby pumps. Omnifuse syringe pumps have been designed to meet the growing needs of the clinical environment. They are suitable for the administration of drugs or other parenteral fluids by or under the supervision of healthcare professionals such as physicians and nurses.

This Instruction Manual describes the features of the standard Omnifuse pump, explains how to use the pump and also describes the pump variants.

The pump has a green Command wheel, LCD surround and keypad.

Omnifuse pump features

The features on the Omnifuse pump can be tailored to suit individual hospital and clinical requirements.

The pump offers a Continuous infusion mode as well as other modes for specialised infusions: Preset Volume, Preset Time, Intermittent and Circadian Rhythm.

The clear messages on the screen guide you through an intuitive programming and running sequence, with further information available from on-screen help explaining how to use a feature or follow a prompt. The Omnifuse pump can be run directly from the AC mains supply or, when necessary, from its own rechargeable backup batteries.

The pump has been designed to be simple to use. Most actions are performed by choosing options from the screen with the Command wheel which you turn and then press. There is a keypad for entry of numerical data, for stopping and starting the pump and for switching it on and off.

Carefully read the entire contents of this manual before using the pump.

Safety and security

Many safety features have been incorporated into the Omnifuse pump:

- Self-test routines run each time the pump is switched on to check it is working correctly.
- Continual checks are also made while the pump is running, especially during an infusion.
- A wide range of audible and visual alarms and warnings, including the LCD screen flashing.
- Each infusion setting must be reviewed and confirmed before an infusion can start.
- Occlusion pressure is monitored as a standard feature. The dry-side pressure level is displayed while the pump is running. The occlusion alarm level setting is adjustable but may be locked to restrict the upper level if required.
- A history of 3000 of the pump's most recent events or actions, stamped with their date and time can be viewed on the pump screen, or transferred to a PC.

User interface

The screen and controls on the pump are easy to use.

- All the information for programming and running infusions is shown on the LCD screen (Liquid Crystal Display).
- Next to the screen is a multi-function Command wheel which you turn and press to select items and enter information.
- The pump has a numeric keypad for entry of data such as bolus dose.
- An Infusion LED flashes to show that the pump is infusing and changes colour if the syringe nears empty. A yellow LED is lit when the pump is connected to the AC mains supply.
- While you program the pump and during an infusion, the screen shows instructions on a message line. Help and information icons lead to more detailed information.
- Status icons show the battery level and occlusion pressure level.
- A special feature lets you program the pump ready for an infusion, then leave it 'asleep' until it is needed.
- The pump totaliser shows totals for the infusion and totals since last reset. A *shift totals* feature is also available.

Administration sets

Graseby medical have a wide range of Flo-Safer administration sets that may be used for specific infusion therapies. These are listed in the *Specification* at the end of this manual.

Pump configuration

The features required for a particular clinical application must be configured on the pump, by a qualified technician, prior to use in the clinical area.

The technician uses Graseby Omnifuse Technician PC software to enable features on the Configuration menu, then transfers this information to the pump. Alternatively, the technician accesses the Configuration menu directly on the pump, using a special password. Once the Omnifuse has been configured for use, it can be programmed with specific infusion settings according to the clinical requirements of the patient and medication.

In addition to the Configuration menu, there is a Technician menu giving access to information about the use of the pump. This menu is normally accessed using PC software but may be accessed directly on the pump, using a special technician's password. For more details about Configuration and Technician menus, read the *Omnifuse Technical User Manual*.

Syringes

On the Omnifuse, the syringe can be loaded as a one-handed operation. Syringes of 2 to 60 ml capacity can be used with the pump.

The syringe type and size are confirmed as the first step in programming an infusion. The syringe types that are available on the pump's menu depend on which have been configured for use.

Bolus features

The Omnifuse pump may be configured to offer one or two methods for administering a bolus: a hands-on and/or a pre-set bolus. The maximum bolus rate can be set up to 800 ml/h although this also depends on how the pump has been set up and the size of syringe in use.

If no bolus is needed for the clinical area, then the bolus feature can be turned off.

Mass Units

A wide range of mass units can be configured for use on the Omnifuse. The pump can be programmed to give a dose of medication which is expressed as a *mass unit*, rather than volume to be infused over a period of time.

Protocols

The pump can offer the clinical user a selection of protocols. Before the pump is used, the hospital/community protocols relevant to the clinical area must be loaded into the pump. This is done using the Omnifuse Drug Protocol Management System.

With Protocols in use, all the user in the clinical area needs to do to set up the pump for a patient is to load the syringe, select the protocol and confirm the parameters displayed on the screen.

Graphics

This feature allows you to view the progress of the infusion in graphic form as an infusion profile.

Power

The Omnifuse is normally run from the AC mains supply. However, it switches automatically to internal rechargeable batteries in the event of a power cut or if it is disconnected from the mains.

On fully charged batteries, the pump may be run for up to 10 hours depending on the infusion rate.

External communications

The pump has an RS232 serial cable connector with PC interface protocol. This is used for configuring the pump from a PC, is used to transfer the pump history to a PC and is also used to control the pump externally.

Omnifuse pump variants

The standard Omnifuse has no In-line occlusion sensing capability or Lockable cover. Three Omnifuse pump variants are available for use in environments where In-line Occlusion Sensing or a Lockable cover are required.

The variants can be specified when the pump is purchased, or the standard pump can be upgraded to one of the variants, if it is returned to Graseby Medical. The three variants are:

- Omnifuse with In-line Occlusion Sensing
- Omnifuse with Lockable Cover
- Omnifuse with In-line Occlusion Sensing and Lockable Cover.

In-line Occlusion Sensing

For clinical applications where occlusion pressure needs to be monitored very precisely, for example in a Neonatal unit, the pump should be ordered with an in-line occlusion sensor.

This is used with a special extension set, which includes a pressure sensing disc, to monitor occlusion pressure in the infusion line while the pump is infusing.

If required, an infusion can be started without the in-line sensor, in which case the pump reverts to the standard dry-side occlusion sensing.

With In-line sensing in use, the Graph hotspot offers a choice of infusion profile graph or a pressure trend graph.

Lockable Cover

The Lockable cover pump variant may be used in clinical areas where it is necessary to prevent the syringe from being removed and to prevent any change to the flow rate for the infusion.

The Lockable cover prevents physical access to the syringe except by an authorised person with a key. The cover can be used with the Omnifuse's Security software enabled, which locks the keypad once the infusion is started. This prevents accidental changes to the infusion rate, or access to the bolus feature.

In-line Occlusion Sensing and Lockable Cover

The Omnifuse is also available as a pump fitted with both the In-line sensor and the Lockable cover.

PC software

With a PC linked to the pump by RS232 cable, you can transfer information between the pump and PC.

Two software packages are available on CD-ROM from Graseby Medical Ltd. for use with the Omnifuse pump.

Omnifuse Drug Protocol Management System

The Protocol support software package is used to design and create Protocols and download them for use on the Omnifuse.

Graseby Omnifuse Technician PC Software

This software package is supplied with the Service Manual. It is required for technical maintenance and servicing an Omnifuse pump, as well as pump configuration. The Technician PC software is used to:

- Perform service tests and calibration on the pump.
- Set parameters and enable or disable options that would otherwise be changed using the pump's own Configuration menu.
- Download and examine the pump's history database.
- Keep a service database for pumps.
- Store configurations and settings from the pump and restore them later.


Omnistack pump stacking system

The Omnistack is a multi-pump stacking system designed for mounting Omnifuse syringe pumps. Up to four Omnifuse pumps can be mounted on an Omnistack, powered from a single mains inlet.

The Omnistack can be clamped to the Graseby wheelbase, or to a suitable pylon. The Omnistack must not be used with a standard IV pole.

For full details on how to use the Omnistack, see the *Omnistack Instruction Manual*.

The part numbers for the Omnistack stacking system, pole and wheelbase are listed in the Accessories section of the Specification at the back of this manual.

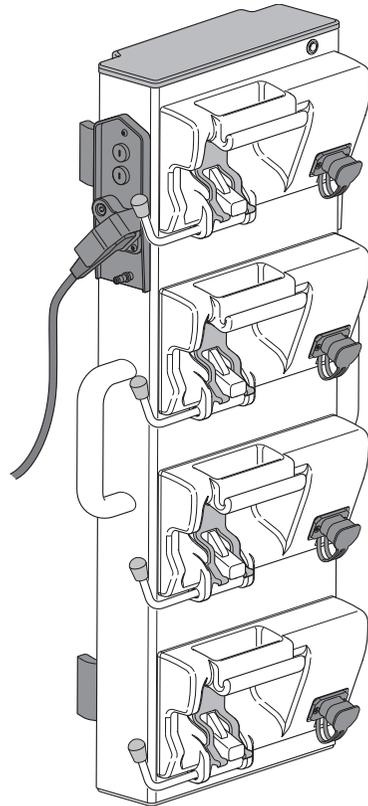


GM19L01001L1

Chapter 2

Basics

	<i>Page</i>
Parts of the Omnifuse pump	1
Switching the pump on and off	4
Omnifuse screen, icons and indicators	8
Loading and unloading a syringe	12
Care and maintenance	18

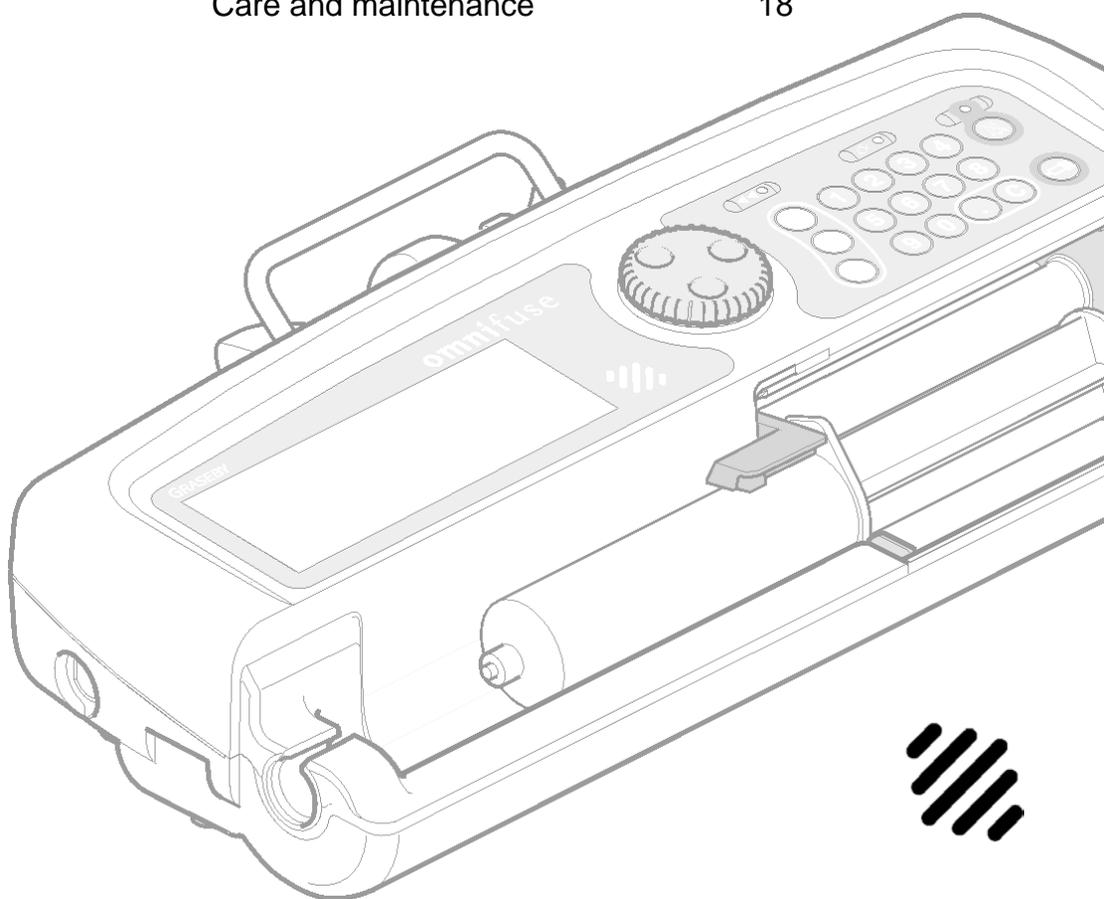


Parts of the Omnifuse pump

This section shows the main external features of the standard Omnifuse pump. It has no In-line occlusion Sensor or Lockable cover.



For information on in-line occlusion sensing on the Omnifuse, see *In-line occlusion sensing, Chapter 6* or for details on the Lockable cover, see *Lockable cover, Chapter 7*.

Rear of the pump

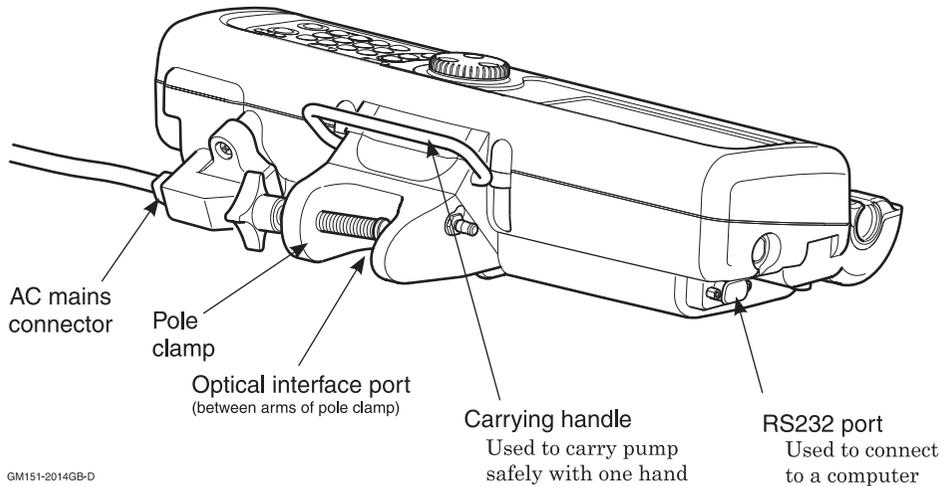

GM151-2014GB-D

Front of the pump

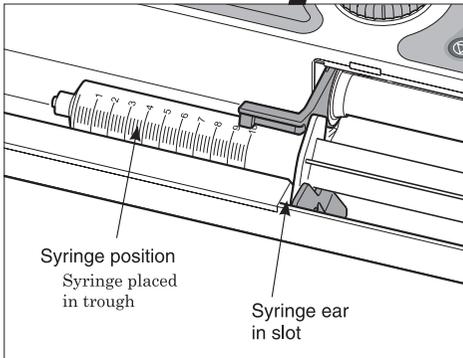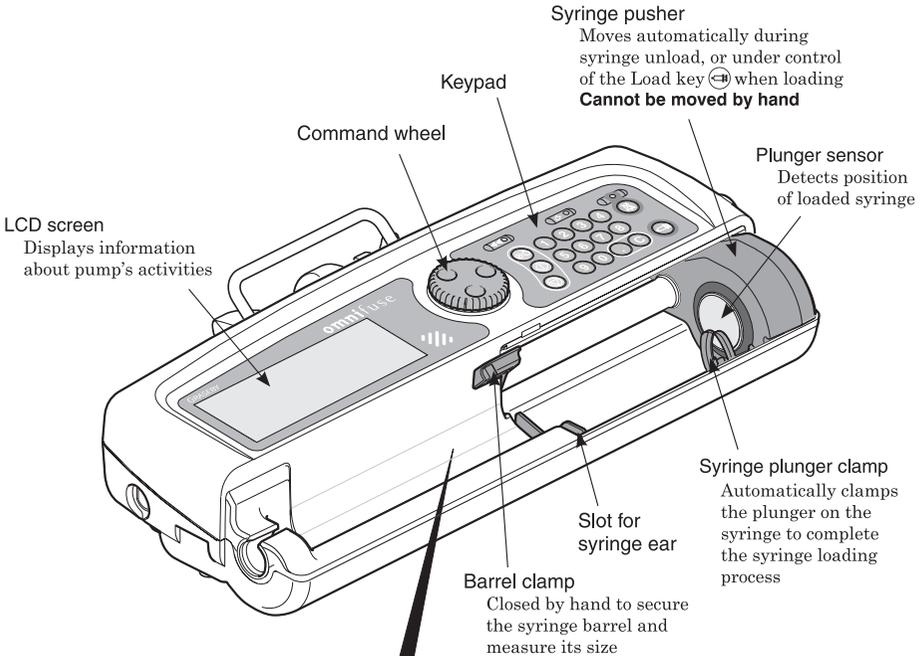

GM151-2012GB-C

This illustration shows the main external features on the front of the Omnifuse pump.

Basics

The Command wheel, keypad and LEDs


Command wheel

You enter commands into the pump with a multi-function Command wheel.

You can turn the wheel freely, either clockwise or counter-clockwise.

If you press the wheel down, you will feel a positive *click*.

- Turn the wheel in either direction to highlight items on the screen such as hotspots
- Press the wheel to select, accept or confirm the parameters or action.

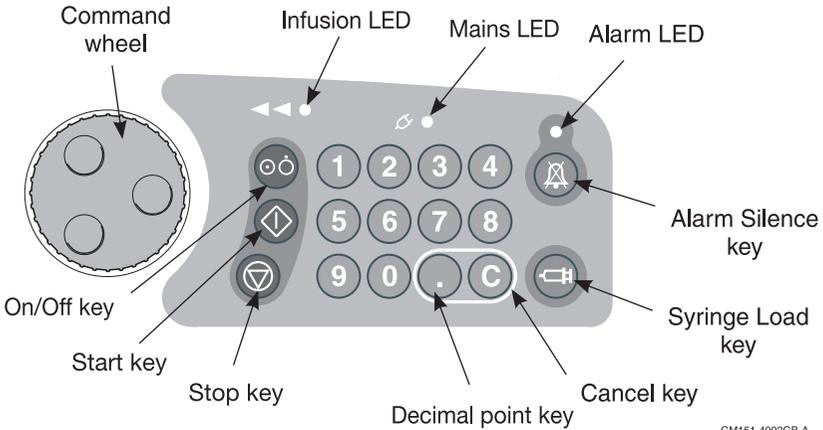

GM151-4002GB-A

Keypad

Next to the Command wheel is the keypad where you will find the:

- ⏻ On/Off key
- ⬆ Start key (Green)
- ⬇ Stop key (Red).

The numeric keypad is for entering numbers. It includes a Cancel key and Decimal point key.

On the right of the keypad is:

- 🔕 Alarm Silence key (Red)
- 📄 Syringe Load key (Blue).

LEDs

The pump has three LEDs:

Infusion LED - This is green, or amber. It is lit or flashes when an infusion is running. For full details, see *Pump indicators*, page 2 - 10.

Mains LED - This is yellow. It is lit when the pump is connected to the AC mains supply and goes out if the pump is operating on batteries.

Alarm LED - This flashes red when the pump is sounding the alarm.

Switching the pump on and off

Before you switch the pump on, visually check for damage to any part of the pump or its connectors.


Plug it in to an AC mains supply if possible. If necessary the pump can be run on its internal batteries see *Using the pump on batteries*, page 2 - 6.

Switching on and power-up tests

To switch on, press On/Off (Ⓞ). The initial screen shows:



The pump continues with its self test:

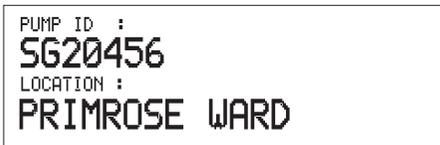


The Pump ID and Location shown on this screen are set up in the Technician menu described in the *Omnifuse Technical User Manual*.

CAUTION: When turning the pump on, if screens similar to those illustrated are not displayed do not use the pump, and send the pump to authorised service personnel.

Service due message

The pump checks its service due date at the end of the power-up tests. If a service is due, you will see a message, for example:

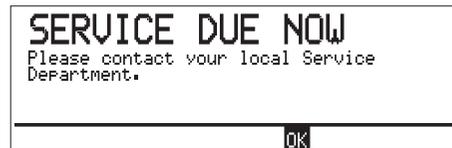

The pump can still be used. To continue with your setup, simply press the Command wheel. If however you wish the pump to be serviced, switch off.

The pump will display the service due message when you next switch it on.

Faulty pump

If the pump discovers a fault during the power-up tests, it displays a message. Do not use the pump, but return it to a qualified service engineer. For more details, see *Troubleshooting, Chapter 8*.

WARNING: Do not use a faulty pump. If the pump detects a fault, a continuous alarm will sound and the screen will display a System Fault message. If this happens, switch the pump off, disconnect it from the mains and take it to a suitably qualified engineer. Incorrect performance of the pump can cause complications resulting in patient injury or death.

Power-up tests completed

When the power-up tests are complete, the screen shows the first of the syringe loading instructions.

For example:



The pump is now ready for you to load the syringe, as described on page 2 - 12.

Switching off

Before you can switch off, the pump must be stopped. If an infusion is running:

- Press and hold the Stop (⏏) key.


The screen displays:



To switch off the pump:

- Press the On/Off key (⏻) and hold it down for two seconds.

As the pump switches off, the screen displays:

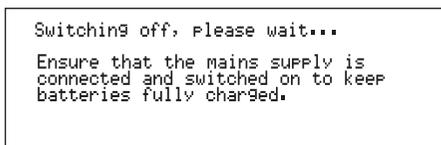

Ensure that the pump is connected to the AC mains to keep the batteries charged and ready for use next time.

Using the pump on batteries

In everyday use, the Omnicuse pump should be connected to a suitable AC mains supply, not used on batteries.

This way, the batteries are kept fully charged, available for use in an emergency.

- The yellow Mains LED is only lit when the pump is connected to the AC mains supply
- The battery icon  is shown on the screen when the pump is running on batteries
- The battery icon is displayed when the pump is connected to the AC mains supply, to show that the batteries are being charged. When the batteries are charged, the icon disappears.

Pump with patient in transit

If a mains supply is unavailable - for example, if the pump is infusing a patient in transit - you can start or continue an infusion with the pump running on batteries.

Mains lead retainer

The mains lead is held in the pump with a retainer. This is so that it cannot be pulled out accidentally.

CAUTION: Failure to use the mains lead retainer means that the pump may be accidentally disconnected from the AC mains supply.

WARNING: Failure to use the mains lead supplied with the pump will compromise the pump's ability to resist fluid ingress, resulting in possible user or patient injury or death.

Mains supply failure

If the AC mains supply fails or if the pump is disconnected from the mains, it automatically switches to its internal batteries.

If the pump is infusing, the infusion continues without interruption. This message is displayed:



- The yellow Mains LED goes out and the battery icon  is shown on the screen when the pump is running on batteries.

The pump will display a warning message approximately 30 minutes before battery power runs out.

When the pump is reconnected to the main supply, or mains power resumes, the pump automatically switches back to the AC mains supply.

Note: When the pump is being used on battery, an independent battery monitor is active. This provides an additional check on the battery condition.

If the pump's alarm sounds but no message is displayed on the screen, the independent battery monitor has detected a problem.

The alarm cannot be silenced and the pump cannot be used. It must be taken to a qualified technician for repair.

Recharging the batteries

While the pump is switched off, you should leave it connected to the AC mains supply so that the batteries are always kept fully charged.

Recharging can take up to 10 hours if the batteries are completely flat.

While the pump is running on the mains supply, the battery icon  is displayed to show that the batteries are being charged. The battery icon disappears once the batteries are fully charged.

Battery life

As a rough guide, an Omnifuse pump can operate on battery power for up to 10 hours at a continuous infusion rate of 5 ml/h with the batteries starting from fully charged. This period includes four syringe load/unload cycles.


To ensure that the batteries perform as specified they must always be fully charged. The batteries must not be left partially charged for extended periods, for example when the pump is in store.

WARNING: Correct management of battery charging is essential to ensure that the pump can operate on batteries for the time specified. Failure to do so may lead to impaired functioning of the pump, resulting in patient injury or death.

Omnifuse screen, icons and indicators

Moving around the screen

To move around on the screen you turn the Command wheel in either direction. The highlighting shows you where you are.

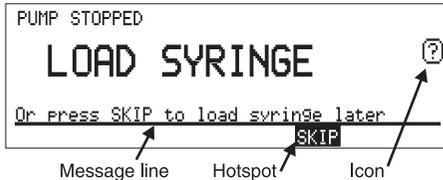

As well as the main command on the screen, you can see:

- The message line
- Hotspots
- Icons.

Message line

This gives you more details about the command on the screen, or shows the status of the infusion, or shows you a warning or alarm message reminder.

Hotspots

A hotspot is an action word on the bottom line of the screen. The hotspots you can see are relevant to whatever the pump is doing. The most common hotspot is **CONFIRM**.


In many cases, the hotspot you need is automatically highlighted on the screen, so you can just press the Command wheel.

If the relevant hotspot is not highlighted, turn the Command wheel to highlight it, then press the Command wheel.

When you press the wheel, the pump will carry out the command. In this example, if you press the Command wheel with **SKIP** highlighted, the pump will skip the syringe load process.

Icons on the pump screen

This example screen illustrates where the Omnifuse screen icons appear:

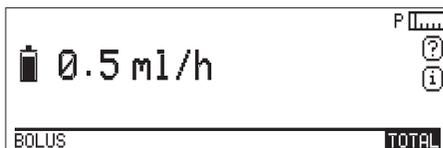

Some icons are displayed on the screen all the time, others only appear in special circumstances.


The table on page 2 - 9 shows all the icons that may appear.

Screen icons

Icon	Name	Description
	Info	Displays settings and parameters for the current infusion.
	Help	Displays a screen of help information relating to the current state of the pump.
	Occlusion Level (Dry-side)	The position of the vertical line shows the level at which the occlusion alarm pressure has been set, from Level 1 on the left to Level 5 on the right. When the pump is running the icon fills from left to right to show the current pressure. You use this icon to set the alarm pressure level, see <i>How to set or change the occlusion alarm level</i> , page 3 - 4.
	Occlusion Level (In-line)	This replaces the dry-side occlusion level icon when the sensor is use on an Omnifuse with In-line occlusion sensing, see <i>In-line occlusion sensing</i> , Chapter 6.
	Battery level	This is displayed as a static icon when the pump is running on batteries, or as a flashing icon when the pump is connected to the AC mains supply and the batteries are charging.
	Hold (Graph)	This appears when the pump is displaying the Graph screen. Select this icon to hold the screen. While the screen is on hold, the icon changes to the return icon shown below.
	Return (Graph)	This is displayed when the pump is holding the Graph screen. Select this icon to return to the running screen, or previous graph screen.
	Zoom (Graph)	This is displayed during long infusions, when the graph display can be zoomed to show more detail for a particular hour.

Pump indicators

There are three LEDs or indicators on the front of the pump, positioned above the numeric keypad.

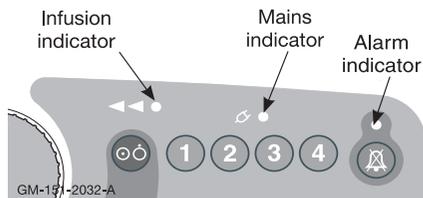

Infusion indicator

The Infusion indicator shows whether or not the pump is infusing, and changes colour when the syringe is nearly empty. See details in the table below.

Infusion LED	Shows ...
<i>Flashing</i>	Medication being delivered.
<i>Steady</i>	Infusion suspended - Stop key (⏏) was pressed.
LED Colour	Shows ...
<i>Green</i>	Syringe is full.
<i>Amber</i>	Syringe is nearly empty.

Mains indicator

The yellow Mains indicator is lit when the pump is connected to the AC mains supply.

Alarm indicator

The red Alarm indicator flashes when the pump sounds the alarm.

Sounds on the Omnifuse

This table describes the sounds that you may hear from the Omnifuse pump and tells you what action is required if you hear them.

For more details on pump warnings and alarms, see *Troubleshooting, Chapter 8*.

Sound	Occurs when...	What you should do
Beep	You press a key or the Command wheel.	If the beep is too soft or too loud, call a technician to change it for you via the Key Beep Volume parameter in the Configuration menu.
Double beep	The pump requires your attention.	Read the message on the screen, then carry out the suggested action.
Alarm	The pump has discovered an operating problem, for example the syringe is empty or there is an occlusion.	Read the message on the screen, then press the Alarm Silence  and carry out the suggested action.
Continuous System Alarm	The pump has detected an internal error.	Read the message on the screen. Pressing Alarm Silence  will have no effect. Switch off the pump and hand it over to a qualified technician for servicing.

Setting the alarm volume

The volume of the alarm on the Omnifuse can be set to one of three levels:

- 1 - Soft
- 2 - Medium
- 3 - Loud.

Provided the pump is not infusing and the cover is open, you can adjust the volume to suit the operating environment.

To adjust the alarm volume:

- Press and hold the Alarm Silence  key and at the same time press the 1, 2 or 3 key on the keypad.

As you press a number, the pump demonstrates the volume of the alarm at your chosen level.

The last number you press sets the volume that the pump will use for the alarm, until set to a different volume.

Loading and unloading a syringe


The pump is designed to make loading a syringe a one-handed operation. There are three components that secure the syringe:

- The syringe ear slot
- The barrel clamp
- The pusher.

The syringe is placed in the trough at the front of the pump. The ear or flange of the syringe should be placed into the syringe ear slot. The barrel clamp is lowered by hand onto the syringe.

The sections on the following pages provide full details on loading/unloading a syringe.

Please note that Warnings relating to syringe loading and unloading are listed on page 2 - 14.

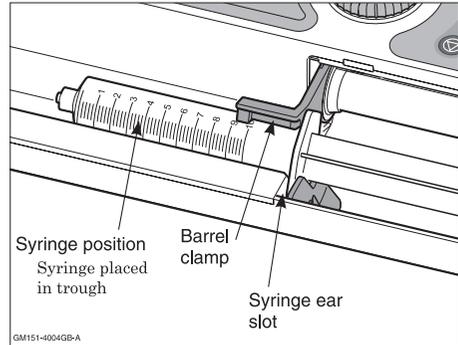

When to load the syringe

When you turn the pump on, you see the screen below, unless a syringe is already loaded:



If a syringe is already loaded you will only be asked to confirm the brand and size.

When to skip loading

If the syringe is not ready, or if someone else is preparing the syringe while you program the pump, you can skip loading. The **SKIP** hotspot allows you to go ahead with programming the pump and load the syringe later on.

Please read *Warnings and Cautions concerning syringes*, on page 2 - 14.

Syringe brand and size

When you load a syringe into the trough and close the barrel clamp the pump displays:

- The name of the previously loaded syringe brand, and
- The size of the syringe it currently detects.

For example:



The pump can sense the *size* of syringe you have loaded, but it cannot sense the *brand*.

Instead it remembers the brand of syringe used last time, then checks the sizes listed for that brand against the size of the syringe sensed.

If the pump finds that the size is not compatible with the brand, it displays a warning on the message line: **Brand does not match detected size**. If you see this message you must select the correct brand. Syringes supported by the Omnifuse are listed in the Specification at the back of the manual.

It is essential that you check the syringe brand displayed by the pump and make sure that it corresponds to the one you have loaded.

Replacing a syringe during an infusion

During an infusion you can change syringes. For example you could change an empty syringe with a full one.

During a Continuous mode infusion

For a pump running a Continuous infusion, you may replace a syringe during the infusion using a different brand, and/or a different size of syringe.


If you need to change the brand, unload the syringe then follow the instructions for *How to change the brand of syringe*, page 2 - 17.

All other types of infusion

The replacement syringe must be the same brand and size as the original syringe brand and size that were

confirmed at the start of the infusion. Suspend the infusion by pressing Stop (⊙), then carry out the instructions for *How to unload a syringe*, page 2 - 16, and *How to load a syringe*, page 2 - 15.

If you try to replace a syringe with a different brand or size, you will see an invalid syringe message:

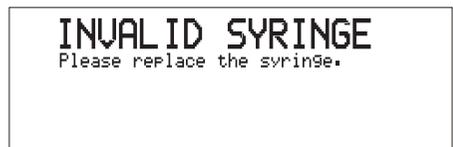

You will not be able to restart the infusion until you load a syringe of the same brand and size as the original.

Please read *Warnings and Cautions concerning syringes*, on page 2 - 14.

Warnings and Cautions concerning syringes

Please read these warnings and cautions before following the instructions in any of these sections:

- *How to load a syringe*, page 2 - 15
- *How to unload a syringe*, page 2 - 16
- *How to change the brand of syringe*, page 2 - 17.

WARNING: Use only the syringes and administration sets listed in Specifications and Standards at the end of this manual. Failure to do so may result in an inaccurate delivery. Graseby Medical does not guarantee performance of the pump if syringes other than those listed are used. Incorrect function or performance of the pump can cause complications resulting in patient injury or death.

WARNING: To avoid syphoning of the syringe contents (free flow), ensure that the syringe is correctly loaded into the pump, that the syringe plunger is properly engaged by the pump's pusher block and that the pump is placed not more than 80cm above the infusion site. Syphoning can result in over-infusion leading to patient injury or death.

WARNING: For safe operation of the pump, the syringe must be correctly loaded. Ensure that the syringe plunger is properly aligned before closing the barrel clamp. Failure to do so may result in inaccurate delivery, resulting in patient injury or death.

CAUTION: Do not attempt to move the pump's pusher by hand. Always use the syringe Load key  to move the pusher. Failure to observe this caution may cause mechanical damage to the pump.

WARNING: The supported syringes are single use only and the administration set should be changed according to the manufacturers instructions. A new syringe and administration set must be used for a new patient. Failure to observe this warning may lead to compromised performance of the pump, resulting in patient injury or death.

WARNING: To avoid over- or under-infusion, always verify that the brand and size of the loaded syringe are the same as the brand and size displayed on the screen before starting an infusion. Failure to do so may result in an inaccurate delivery of medication, resulting in patient injury or death.

WARNING: To avoid syphoning of the syringe contents (free flow), ensure that the patient line is clamped before loading or unloading the syringe. Syphoning can result in over-infusion leading to patient injury or death.

WARNING: Ensure that your fingers are not in the path of the pusher during syringe loading or unloading. Failure to do so may result in user injury.

How to load a syringe



1. Ensure that the line to the patient is clamped, or disconnected. Switch on the pump if it is switched off. Lift up the barrel clamp if it is closed - you will see this message
2. Place the syringe in the trough. The syringe ear must fit into the slot.
3. Lower the barrel clamp by hand.
The screen shows the size of syringe the pump has detected and the brand of syringe used last time.
4. Press the Command wheel to **CONFIRM** that the syringe brand and size are correct. To change brand, see *How to change the brand of syringe*, page 2 - 17.
5. Press and hold the blue Load key .
6. Keep positive pressure on the Load key  while the screen shows this message.
7. Release the Load key  when the the pump beeps and displays this message. The pump will display the next screen in the programming sequence.

Please read *Warnings and Cautions concerning syringes*, on page 2 - 14.

How to unload a syringe

If the pump is switched off with a syringe still in place, you do not need to turn it on: lift the barrel clamp and gently remove the syringe from the trough.

If the pump is switched on, the infusion must be suspended, or the pump must be stopped before you try to unload a syringe.

1. Press Stop . If you press the key briefly, the infusion is suspended, or if you press and hold Stop  the pump displays: **PUMP STOPPED**.
2. If infusate is still present in the syringe, clamp the patient line before proceeding.




3. Lift up the barrel clamp by hand.

The pump recognises that you have lifted the clamp and asks you to confirm that you wish to remove the syringe.

4. Confirm **YES** by pressing the Command wheel.

The pusher will automatically release the syringe plunger and retract.

5. When you see this message, remove the syringe from the trough.

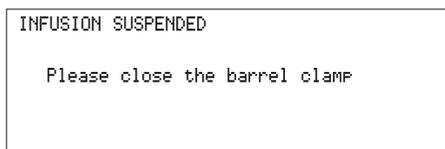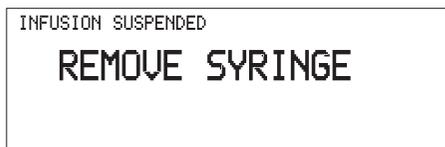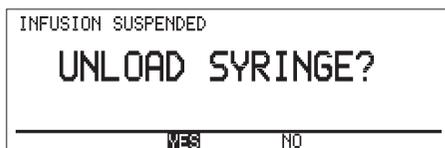

To discontinue unloading:


If you did not intend to unload the syringe:

1. At Step 4 above, turn the Command wheel to select **NO**. The pump asks you to close the barrel clamp.
2. Close the barrel clamp. Press Start  to resume the infusion.

Please read *Warnings and Cautions concerning syringes*, on page 2 - 14.

How to change the brand of syringe

If the syringe brand and size displayed are not correct when you load the syringe, follow the steps shown here:

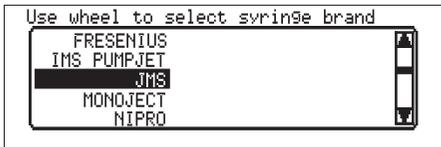

1. Turn the Command wheel to highlight the syringe brand, then press the Command wheel to change the selection.
2. The screen shows a menu of syringe brands. Turn the wheel to highlight the correct brand, then press the Command wheel to select it.
3. When the syringe brand and size are correct, press the Command wheel to **CONFIRM**.
4. Press and hold the blue Load key .
5. Keep positive pressure on the key while the screen shows this message. When the pusher is properly engaged with the syringe plunger, the pump beeps.
6. When the pump displays this message, release the Load key . The pump will display the next screen in the programming sequence.

Please read *Warnings and Cautions concerning syringes*, on page 2 - 14.


Care and maintenance

This section explains the everyday care of the Omnifuse pump. It does not cover technical pump maintenance, which is

described in the *Omnifuse Technical User Manual*, and the *Omnifuse Service Manual*.

Carrying the pump

The Omnifuse pump has a carrying handle, which you should always use to carry the pump safely.

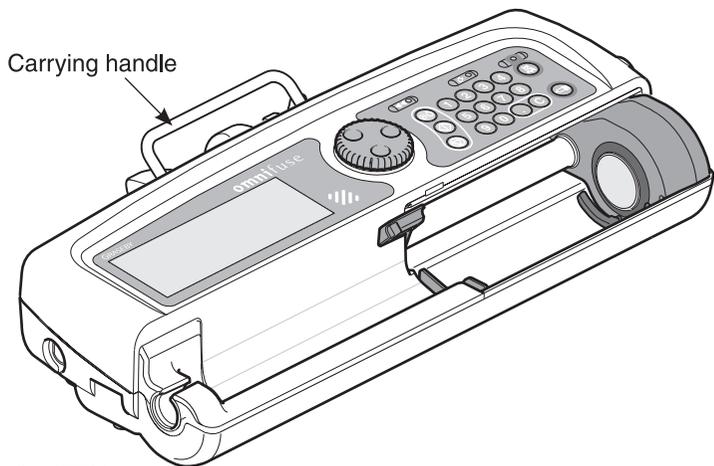

CAUTION: Never carry the pump except by the handle. Failure to do so may result in damage to the case, or you may drop the pump and cause it internal damage.

GM151-2003GB-D

Positioning the pump for use

When you position the pump in a location suitable for the patient, make sure it is:


- Clamped to an IV Infusion pole
or
- On a stable horizontal surface such as a table where it cannot slide or slip. The surface should be at least 500 mm x 250 mm to accommodate the pump
or

- Fitted into an Omnistack pump stacking system, which can hold up to four Omnifuse pumps.

The Omnistack must be used with the Graseby wheelbase, or a suitable pylon. See the *Omnistack Instruction Manual* for full details.

Using the pole clamp


The Omnifuse pump is fitted with a pole clamp designed to hold the pump at a 45° angle against an IV infusion pole.

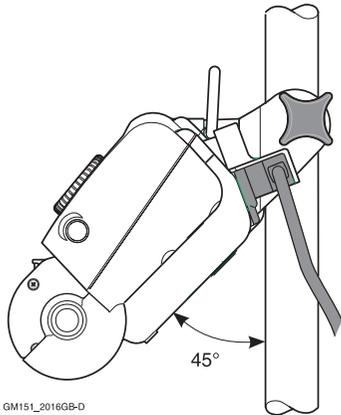

This angle makes it easy for you to read the screen and also gives you good access to the keypad and Command wheel. To clamp the pump to an IV pole:

1. Open the clamp screw so that the clamp is open wide enough to slip onto the pole.
2. Hold the pump with both hands and position it with the pole between the jaws of the clamp.
3. Supporting the pump underneath with one hand, tighten the clamp screw with the other hand until the pump is securely fixed to the pole.

The rear of the pump looks like this with the pole clamp correctly fixed to a pole.

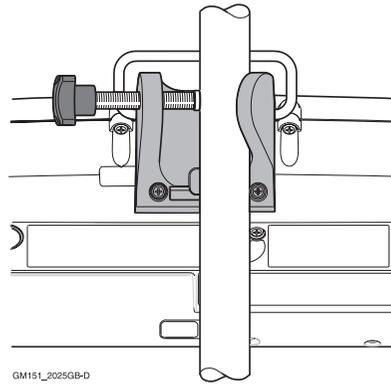

WARNING: To avoid the pump becoming detached from an IV pole always make sure that the pump is securely fixed to the pole. Always check the security and stability of the assembly with the pump mounted. If no IV pole is used make sure that the pump is completely stable on a horizontal surface. Failure to observe this warning may cause damage to the Omnifuse pump and harm the operator or the patient. As a result, the operator or patient may suffer direct injury, or the Omnifuse pump may fail to operate correctly, leading to patient injury or death.

Cleaning and care of the pump

Cleaning the pump

The Omnifuse pump is designed to be water-resistant against accidental spillages, but not waterproof. You should therefore clean the casing and outer surfaces of the pump using a damp cloth, or if necessary a cloth dampened with a mild solution of washing-up liquid.

Clean the pump as follows:

1. For Omnifuse pumps with a Lockable cover, remove the cover as described in *Removing and replacing the cover*, page 7 - 5, then clean the cover.
2. Wipe over the exterior surfaces of the pump, paying particular attention to the barrel clamp and the syringe ear slot. To clean beneath the syringe pusher block, switch on the pump and load an empty syringe so that the block moves left.
3. When cleaning is complete, disinfect the pump using a suitable disinfectant solution and remove any disinfectant residue by wiping with a clean damp cloth.
4. After cleaning and disinfection, remove the syringe and replace the Lockable cover.

WARNING: Following a significant liquid spill onto the pump, it should be wiped dry and inspected by service personnel before being returned to service. Failure to do so may result in compromised functioning of the pump, leading to patient or user injury or death.

CAUTION: Do not use cleaning and disinfecting agents other than the approved ones specified here.

Spray cleaners

If using a spray cleaner, do not spray the pump, but dampen a cloth with the cleaner and use this to wipe the pump.

Disinfectants

Disinfect the casing of the pump either with a cloth dampened in a solution of sodium hypochlorite (0.2%), or with alcohol wipes intended for disinfecting equipment.

A suitable disinfectant solution can be made by diluting sodium hypochlorite with water to give a solution of 0.1% available chlorine. Preferably use a freshly made solution, and do not use one which is more than 24 hours old.

Maintenance

Other than cleaning and disinfection, the Omnifuse pump requires no maintenance to be carried out at ward-level. In particular it does not require lubrication - and indeed any lubricant containing an organic solvent may damage the plastic of the pump casing.


If you consider that the pump needs further attention after you have cleaned and disinfected it, return it to a qualified engineer for servicing.

CAUTION: The pump must not be immersed in any liquids or exposed to strong organic solvents. Wipe off spills immediately, and do not allow fluid or residues to remain on the pump. Additionally, the pump is not designed to be autoclaved, steam-sterilised, EtO-sterilised or subjected to temperatures in excess of 45° C (113° F). Failure to observe this caution may cause serious damage to the pump.

Chapter 3



Programming the pump

	<i>Page</i>
Continuous infusion in ml per hour	1
Occlusion pressure and alarms	3
The Totaliser	6
Purging the line	10
Sleep mode	11
Omnifuse bolus feature	13
Graphics	17

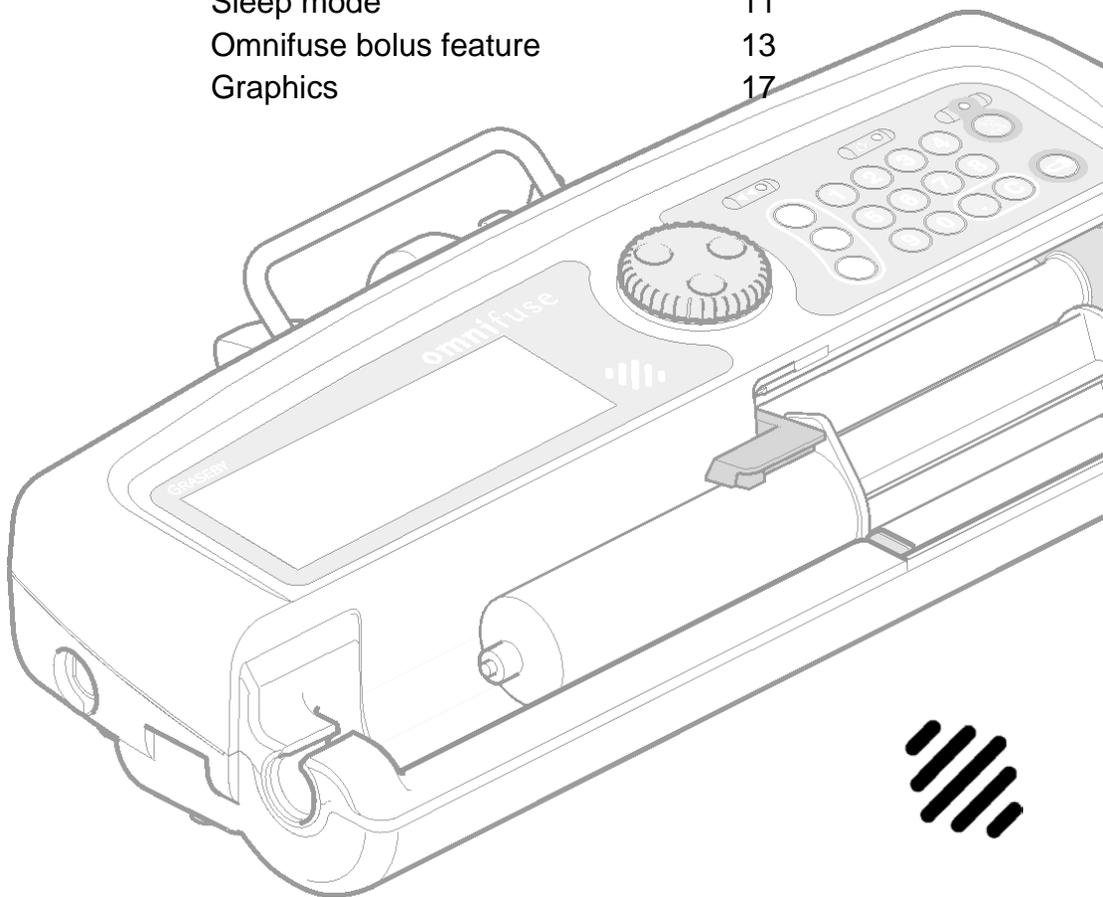

Continuous infusion in ml per hour

This is an example of the screen you will see once you have loaded the syringe on a pump which is used for continuous infusions only:

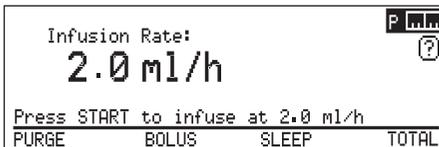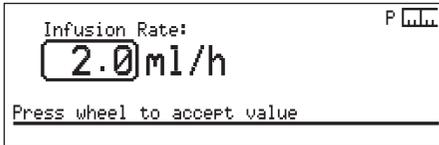

1. Check the infusion rate displayed on the screen.
2. If you need to adjust it, enter a new rate from the keypad or turn the Command wheel.
3. The infusion is ready to start at this point, but you should first:

When the rate is correct, press the Command wheel to accept the value.

- Check or change the Occlusion Alarm Level to ensure that it is appropriate for the infusion, see page 3 - 4.
- Select **TOTAL**, to check or reset the cumulative total, see page 3 - 6.

You could also carry out these tasks if necessary:


- **PURGE** the line see page 3 - 10.
- Use the **SLEEP** feature see page 3 - 11.
- Give a **BOLUS** see page 3 - 15.

4. When you are ready to start the infusion, press the green Start key .
5. The running screen shows the infusion rate, and the total volume infused since the last reset. The Infusion LED flashes to show that the pump is running.

Changing the rate while running

Once a Continuous infusion is running, you can alter the infusion rate if necessary.

The *Rate change while running* feature can be turned off by a technician. If so, you will not be able to highlight the Infusion Rate.

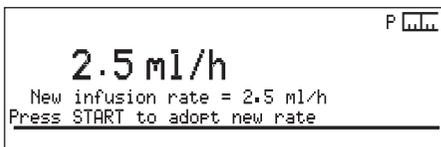

1. Highlight the current rate on the running screen and press the Command wheel.
2. Turn the Command wheel, or use the numeric keypad to enter the new rate.
Press the Command wheel to accept.
3. Press Start  to continue the infusion at the new rate.
If you do not press Start , the infusion continues at the previous rate, and the pump displays a warning message to remind you that the Rate change was not completed.

Ending a Continuous infusion

To end a Continuous infusion:

- Press and hold Stop  - you must keep it pressed for at least two seconds.

If you only press the Stop key briefly, the infusion will be suspended, and can be restarted with the Start key .

Occlusion pressure and alarms

This explains how to monitor the occlusion pressure using the standard Omnifuse dry-side pressure display. Dry-side infusion pressure is detected by a sensor inside the pusher block as it presses against the syringe plunger.

The actual pressure in the infusion line depends on a number of factors such as the viscosity of the liquid, the stiction of the syringe plunger, the bore of the line and the height of the pump relative to the patient. These variables mean that the pressure of liquid in the line can be measured only approximately at the syringe plunger.

For more precise pressure monitoring, an Omnifuse with in-line pressure sensing is available, see *In-line occlusion sensing, Chapter 6*.

Omnifuse pressure levels


You set the dry-side occlusion alarm in *levels*. There may be up to five levels ranging from Level 1 (low) to Level 5 (high). The lower you set the alarm level, the faster the response time will be to a potential problem with the infusion line. For example, if Level 1 is selected, then as soon as pressure in the line begins to build up, the alarm will sound.

Note: If necessary, the pump will raise the occlusion alarm level temporarily during a bolus, to avoid a nuisance occlusion alarm. You may notice the change in the pressure icon.

At the end of the bolus, the alarm level automatically returns to the selected level.

The occlusion alarm is slow to sound if Level 5 has been selected, so this could be used if the patient is under close supervision, so that there is less need to rely on the alarm.


The Occlusion Levels menu shows the levels available, and the approximate pressure for the alarm to sound for each level.

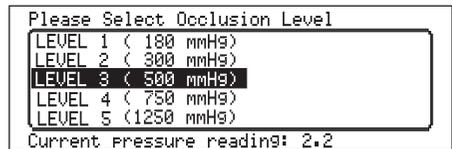

Occlusion alarm lock level

In some clinical areas it may not be safe for the occlusion alarm to be set to a high level. In this case, a technician can configure the pump so that the menu only shows the levels that are acceptable.

For example, the screen looks like this if the Occlusion Alarm Lock level has been set to 3:

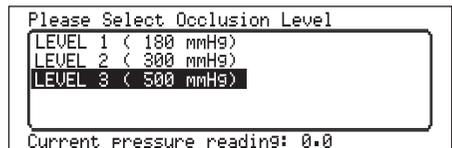


For full details of how to lock the occlusion alarm level, see the *Omnifuse Technical User Manual*. To find out how to change the occlusion level and monitor the alarm pressure during an infusion, see page 3 - 4.


How to set or change the occlusion alarm level

You should always check the occlusion alarm level, to ensure that it is appropriate, before you start an infusion. You can also check or change the level while the infusion is running.


Follow these instructions to check and change the occlusion alarm level:

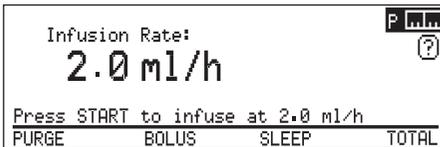


1. Highlight the pressure icon:  and press the Command wheel.

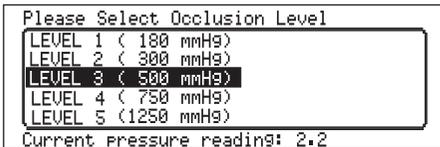

2. The Occlusion Level menu is displayed, showing the currently selected level.

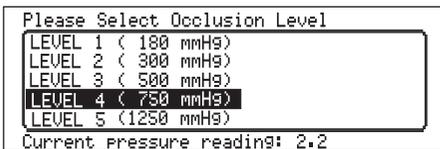

3. Turn the Command wheel to highlight the new level at which you want the occlusion level alarm to sound.

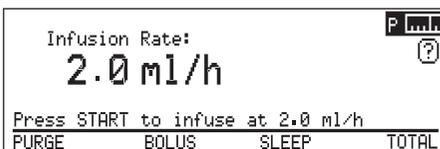

4. Press the Command wheel to set the level and exit from the Occlusion Level screen.

The pressure icon now shows the new level.

WARNING: If an occlusion alarm occurs, immediately clamp the line to eliminate the possibility of a bolus being delivered to the patient. Then inspect the fluid pathway for kinks, clogged catheter, etc. in order to remove the occlusion prior to restarting the infusion. An unintentional bolus of medication can result in patient injury or death.

WARNING: The occlusion alarm level must be checked before starting an infusion to ensure that it is appropriate for that infusion. Failure to do so may result in an unacceptably slow time to occlusion alarm, resulting in patient injury or death.

Monitoring the dry-side pressure

The icon  shows the selected occlusion alarm level, and a visual guide to the current infusion pressure in the line. The line here:  shows that it is set at Level 3.

While an infusion is running the pressure icon fills from left to right to indicate a rise in pressure.

In the icons shown in the example below, the alarm is set at Level 5:



You monitor infusion pressure by seeing where the alarm level has been set and how close the actual pressure is to it.

The pressure level will normally fluctuate as the infusion progresses. If you notice that it is steadily moving closer to the alarm level, you should check the patient and line immediately to find out why.

Note: The standard dry-side occlusion sensing offered on the Omnifuse is adequate for most purposes.

If you require more sensitive monitoring of infusion pressure, with faster warning of an occlusion, you should use one of the pump variants which have an In-line Occlusion Sensor, see *In-line occlusion sensing, Chapter 6*

The Totaliser

The Omnifuse pump has a Totaliser feature so that you can monitor the volume of fluid infused by the pump. The pump can show three separate types of total:

- Cumulative total
- Hourly total
- Shift total.

The totals include any bolus given, and the bolus volume is also displayed separately.

Note: Cumulative and Hourly totals are always available. Shift totals are shown in the example screens here, but may not be in use on your pump.

About the cumulative total

The cumulative total is displayed on the message line when the pump is running, and can also be viewed from the **TOTAL** hotspot.

```
Total delivered since reset =
  20.0 ml (Bolus :   2.0 ml)

Totals reset at : 13:29 03/02/03
THIS HOUR THIS SHIFT  RESET  HOLD
```

The cumulative total needs to be reset manually when required - it is not reset automatically. It may be used for monitoring fluid balance for a patient.

You can reset the cumulative total from the **TOTAL** hotspot:





- Before you start the infusion, for example if the programmed infusion is for a new patient
- Once the infusion has started. In this case, you must press Stop (⊙) to suspend the infusion before using the **TOTAL** hotspot.

Very long infusions

If you run an infusion that continues beyond 120 hours (or five days), then the oldest hour-by-hour data will be overwritten. This means that you can view no more than the last 120 hours. However, this does not affect the cumulative infusion total which will continue to increase until it reaches the limits of 9999 ml for volume and 99999.99 mg for drug mass.

How to view and reset the cumulative total

This example shows how to view and reset the cumulative total. You can do this before starting the infusion, or if the infusion is suspended.

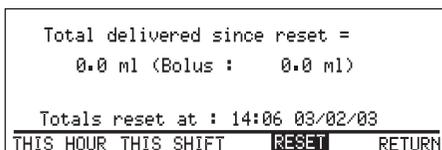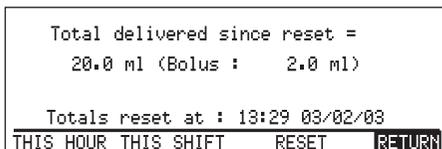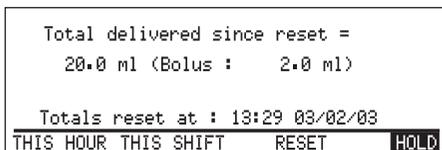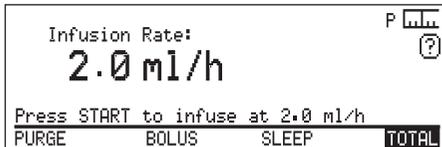

1. On the Infusion Start screen, highlight the **TOTAL** hotspot, then press the Command wheel.

The Total delivered since the last reset is displayed. The totals will be displayed on the screen for five seconds.

2. To view the display for longer, highlight the **HOLD** hotspot and press the Command wheel.

On the total screen, you can see the:

- Cumulative total. This will also show total drug mass if the dose and/or rate was set using mass units.
- Bolus total since the last reset.
- Time and date of the last reset.
- Hotspots leading to the other types of total: **THIS HOUR** and **THIS SHIFT**, see page 3 - 8.

3. To reset the total, highlight **RESET** then press the Command wheel.

4. Select **RETURN** when you have finished, to go to the Infusion Start screen.

Hourly total and Shift total

Hourly totals and Shift totals work in a similar way. The shifts must be defined and programmed into the pump for a specific clinical environment and may be switched off or changed with the Configuration menu, or Graseby Medical Technician's PC software, see the *Omnifuse Technical User Manual*.

The Hourly and Shift totals are reset automatically at the start of each infusion.

They are accessed from the **THIS HOUR** and **THIS SHIFT** hotspots on the Cumulative Total screen.

Note: The Hourly totals are always available, but Shift Totals can be turned off on a pump. In this case, the **SHIFT TOTAL** hotspot does not appear on the Cumulative total screen. See the *Omnifuse Technical User Manual*.

```
Total delivered since reset =
  20.0 ml (Bolus :   2.0 ml)

Totals reset at : 13:29 03/02/03
THIS HOUR THIS SHIFT RESET HOLD
```

```
Total delivered this hour =
  10.0 ml (Bolus :   1.0 ml)

Since : 13:00 03/02/03
PREVIOUS TOTAL HOLD
```

```
Total delivered this shift =
  10.0 ml (Bolus :   1.0 ml)

Since : 18:30 03/02/03
PREVIOUS TOTAL HOLD
```

1. On the Cumulative total screen, turn the Command wheel to highlight **THIS HOUR** or **THIS SHIFT**, then press the Command wheel.
2. **THIS HOUR** displays a screen showing the total infused in the current hour.

THIS SHIFT shows the total delivered since the start of this shift.


For longer infusions, you will see a **PREVIOUS** hotspot. To view the total for longer, select **HOLD**.

You can add together **THIS HOUR** and **PREVIOUS** figures to find out the total volume infused since the start of the infusion, see page 3 - 9.

3. To return to the Cumulative total, highlight **TOTAL** and press the Command wheel. To go back to the running screen, leave the keypad untouched for a minute or press **RETURN**.

Understanding the hourly total

To understand **THIS HOUR** and **PREVIOUS** totals, look at this diagram:



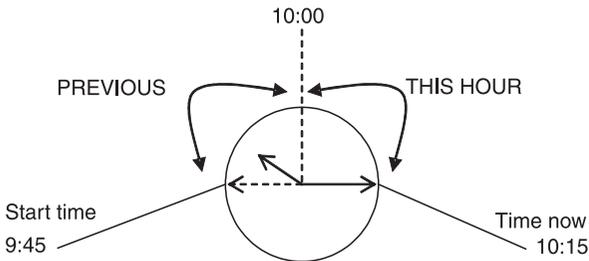


If you pressed Start  at 9:45am to commence an infusion, and checked the totals at 10:15:

- **THIS HOUR** would show you the amount delivered since 10:00. The **PREVIOUS** hotspot would appear on the screen.
- **PREVIOUS** would show the volume if used between 9:45 and 10:00.
- **TOTAL** would return you to the cumulative total screen.

Purging the line

The purging of the line may have been completed before loading the syringe. However, purging with the **PURGE** facility on the Omnifuse pump ensures that the take-up time is reduced to the shortest possible time.

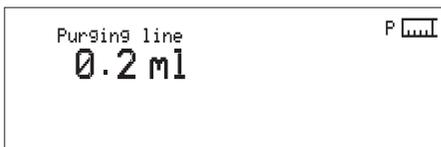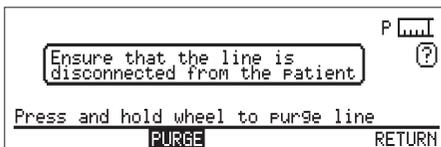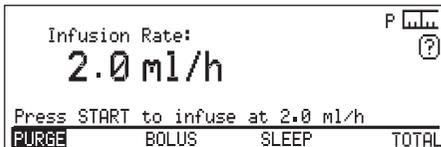

WARNING: To avoid patient embolism, ensure that the patient tubing is purged of all air bubbles before administering any medication. The pump provides a purge facility to assist with this process. The presence of air within the medication can result in complications leading to patient injury or death.

The purge rate is set by a technician, in the range 50-800 ml/h. Details of how to set the purge rate are in the *Omnifuse Technical User Manual*.

1. Highlight the **PURGE** hotspot and press the Command wheel.
2. The pump displays a warning to ensure that the line is *not* connected to the patient.
3. The **PURGE** hotspot is highlighted. Press and hold the Command wheel to purge the line.
4. The pump will beep as it purges and the volume that has been purged through the line will be displayed. You may release the Command wheel at any time and the purge will stop. The pump will automatically stop purging when it reaches 2 ml on the screen.

WARNING: The volume of fluid contained in the connecting tubing is a residual amount and will not be infused. Allowance must be made for this extra volume of fluid when initially filling the syringe and purging the system. Under-delivery of medication can cause complications resulting in patient injury or death.

WARNING: To avoid over-infusion, do not purge the infusion line when the administration set is connected to the patient. Over-infusion of medication can result in patient injury or death.

Sleep mode


You may want to program the pump ready for an infusion, but connect it to the patient and start the infusion later.

For example, you may program a pump while the patient is in theatre, so it is ready to use when the patient arrives in the recovery area.

In this case, you should program the pump, then use the Omnifuse's standby mode, called **SLEEP** mode. In this mode, the pump can be left switched on and ready to start, but does not give any **Not Infusing** warnings.

How to use SLEEP mode

To put the pump into sleep mode:

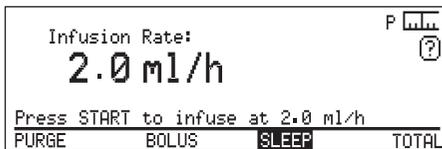

1. Program the infusion and confirm all the infusion parameters. The Infusion Start screen will display the **SLEEP** hotspot.

Turn the Command wheel to highlight **SLEEP**

2. Press the wheel. The pump displays this message.

The pump will stay **ASLEEP** until you press the Command wheel to wake it up.

While the pump is asleep

While it is asleep, the pump does not sound the Not Infusing alarm.



All the other alarms, for example, the low battery alarm, or syringe tampering alarm will still sound to warn you about conditions on the pump.

You can lift the barrel clamp to load or change the syringe when the pump is asleep, without waking it up.

The screen shows the usual syringe loading messages then displays the **ASLEEP** screen again once the loading is completed.

Waking the pump up

When you are ready to connect it to the patient and start the infusion, you need to wake the pump up.

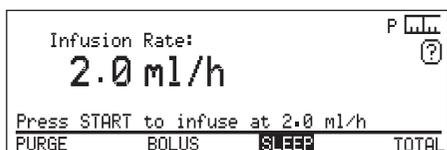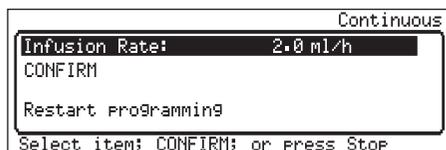

1. The **WAKE ME** hotspot is already highlighted. Press the Command wheel to wake the pump and exit from Sleep mode.
2. The pump takes you directly to the Review screen for the infusion that was programmed earlier.
If you need to change settings at this point, highlight a setting and press the Command wheel.
3. When the settings are correct, turn the Command wheel to highlight **CONFIRM** and press the wheel to display the Infusion start screen.

Omnifuse bolus feature

The **BOLUS** hotspot appears whenever you can deliver a bolus:

- Before starting the infusion, once you have confirmed all the infusion parameters (Continuous mode only)

- During the infusion or when the infusion is suspended.

If the **BOLUS** hotspot does not appear at all, then the feature has been turned off by a technician.

Types of bolus delivery methods

You may have just one, or a choice of two bolus delivery methods, depending on how the technician has configured the pump. The possible bolus methods are:

- Hands-on
- Pre-set.

Hands-on bolus

You keep the Command wheel pressed down until the pump has delivered the bolus dose required.

Pre-set bolus

You enter the bolus amount, then press the Command wheel to start the bolus. Once started, you do not need to hold the Command wheel pressed down.

When it has delivered the dose, the pump automatically ends the bolus.

Default rate for bolus

The default rate for a bolus can be up to 800 ml/h, although it will be lower if the pump is loaded with a small size of syringe. The default rate is set by a technician.

You can change the rate before delivering the bolus.

Bolus beep

The pump will beep once every second while you deliver a bolus. This feature can be turned off by a technician.

WARNING: Where the bolus function is not required, it should be disabled to prevent inappropriate use. Failure to observe this precaution may result in patient injury or death.

Quitting the Bolus screen

If you select **BOLUS** but decide not to go ahead, do not touch the Command wheel or keypad for 10 seconds. The pump will then return to the previous screen.

Bolus and TOTAL

The Totaliser shows the bolus volume as a separate total. To view this, select **TOTAL**:

```
Total delivered since reset =
  20.0 ml (Bolus :    2.0 ml)


Totals reset at : 13:29 03/02/03
THIS HOUR THIS SHIFT  RESET  HOLD
```

The total delivered since reset includes the fluid given as a bolus. In this example:


- The total fluid delivered since the last reset is 20 ml
- The amount delivered as a bolus that is included in this total is 2 ml.

See *The Totaliser*, page 3 - 6 for further details.

How to give a hands-on bolus

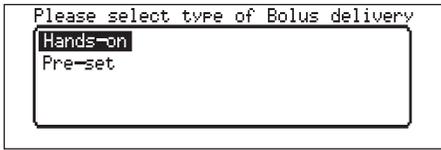


1. Highlight the **BOLUS** hotspot and press the Command wheel.
If this menu appears, highlight **Hands-on** and press the Command wheel.

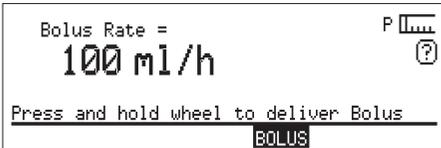


2. The screen will show the bolus rate previously programmed. If this is applicable for the patient, press and hold down the Command wheel to deliver the bolus.

If not, see *Changing the bolus rate*, below.

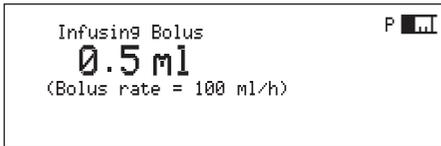


3. Watch the screen to see the amount delivered so far. When the dose is complete, release the Command wheel.

Changing the bolus rate:

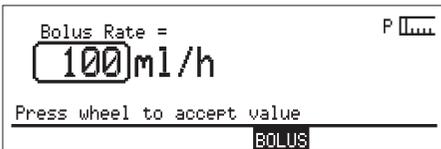


1. If the bolus rate displayed is not applicable, turn the Command wheel to highlight the bolus rate, then press the Command wheel.

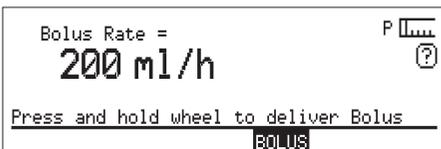

2. Change the bolus rate using either the Command wheel or numeric keypad.


Press the Command wheel to accept the new rate.

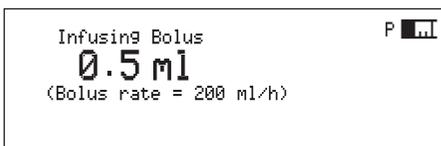

3. Press and hold down the Command wheel to deliver the bolus. When the dose is complete, release the Command wheel.

How to give a pre-set bolus

Pre-set may be available as an additional bolus delivery method.

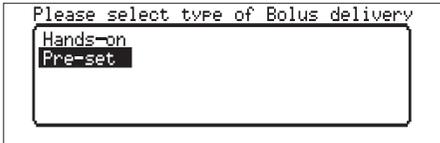

1. Highlight the **BOLUS** hotspot and press the Command wheel.

2. On the Bolus delivery menu, highlight Pre-set and press the Command wheel.

3. Enter the dose required using the Command wheel or the numeric keypad. The pump shows the time required to deliver the bolus dose.

4. Press the Command wheel to accept the dose.

If necessary you can change the rate by adjusting the rate or time, or you may change the dose:

- Turn the Command wheel to highlight the field, press the Command wheel and then enter the new value. Press the Command wheel to accept.

5. If the settings are applicable for the patient, start the bolus:

- Press and hold the Command wheel for one second, then release it - the bolus is delivered automatically.

6. If you need to stop before the whole dose has been delivered, press the Command wheel.



Graphics

Using the Omnifuse Graphics feature, you can view the progress of an infusion as an Infusion Profile graph.

The graph is shown with annotations for alarms, bolus and syringe changes.


For infusions lasting more than an hour, a zoom tool is provided for closer examination.

A text screen summarises the information displayed on the graph.

Note: You can display a graph showing the In-line pressure trend from the Graph hotspot on the Omnifuse variant with In-line Occlusion Sensing, see *In-line occlusion sensing*, Chapter 6.


Accessing the Graphics screen

The **GRAPH** hotspot is displayed during an infusion.

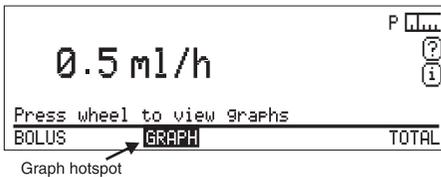

Highlight the **GRAPH** hotspot and press the Command wheel.

The screen shows an **Infusion profile** graph looking something like this:

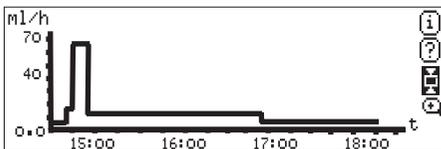

For details on how to interpret the information on this screen, see *Viewing infusion data as a graph*, page 3 - 18.

Viewing infusion data as a graph

The Graphics feature can display the previous 24 hours of infusion data. If the infusion lasts longer than 24 hours, the oldest data is lost.

As the infusion proceeds the displayed time period changes to show either the

last hour, the last four hours, the last 10 hours or the last 24 hours.

If the infusion lasts for more than an hour, you can examine the profile for each hour in greater detail by zooming in on it, as described below.

Graph icons

The graph screens show icons in the top right-hand corner:

Icon	Name	Description
	Info	Displays settings and parameters for the current infusion.
	Help	Displays a screen of help information relating to the current state of the pump.
	Hold	Select to hold the screen display for one minute. While the screen is on hold, the icon changes to the Return icon below.
	Return	This is displayed when the pump is holding the Graph screen. Select this icon to return to the running screen, or previous graph screen. You can also use it to return to the full profile when you have zoomed in on a particular hour of a long infusion.
	Zoom	This is displayed during infusions longer than one hour, when the graph display can be zoomed to show more detail for a particular hour. Use the Return icon to go back to the full profile.

Graph annotations

Letters are used on a graph to indicate special events at the time they occurred:






- **A** means an alarm
- **B** means a bolus
- **S** means a syringe change

If an alarm and a bolus are close together the pump will try to display both letters. However, with long infusions this may

not always be possible. In such cases an alarm annotation will always be shown in preference to a bolus annotation. You will normally be able to recognise a bolus anyway since the rate will change relative to the rest of the profile. Use the Zoom icon to examine that part of the graph in more detail.

How to view an Infusion Profile graph and text

This section explains how to view an Infusion Profile graph and see the text information available from the graph.

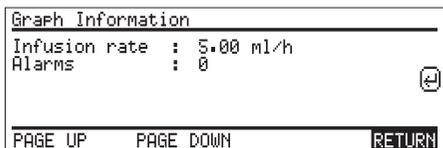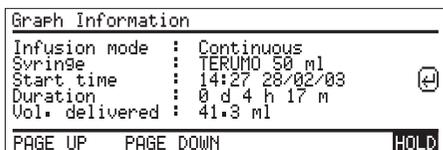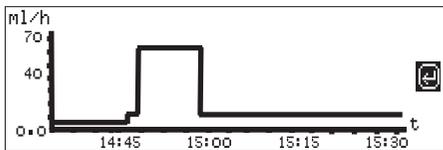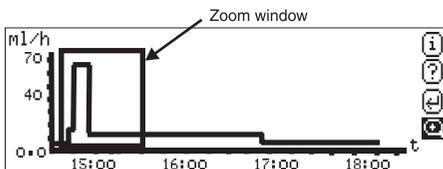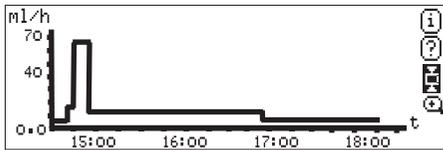

1. Select the **GRAPH** hotspot to display the infusion profile graph, see page 3 - 17.

Select the Hold icon  to keep the graph on the screen for a minute.


2. For a long infusion, select the Zoom icon  to display the zoom window.
3. When the Zoom window appears, turn the Command wheel to move the window left or right along the graph until it is positioned over the hour that you want to examine.
4. Press the Command wheel to display the selected hour in greater detail.
5. To leave the graph, select the Return icon .
6. To view the text screen associated with the graph select the Info icon . This displays a page of text relating to the infusion.
7. Press the Command wheel to **HOLD** the screen for one minute so that you can read the first page of details.
8. Select **PAGE DOWN** to see the second and subsequent pages. **PAGE UP** takes you to previous pages.
9. When you have finished select Return icon  to go back to the graph screen, or **RETURN** to leave the Graphics display altogether.

Programming the pump

Chapter 4

Specialised infusions

	<i>Page</i>
Mass units and special infusion modes	1
Programming using mass units	2
How to change the infusion mode	7
The Review screen	8
Preset Volume infusion	9
Preset Time infusion	12
Intermittent infusion	15
Circadian Rhythm infusion	19

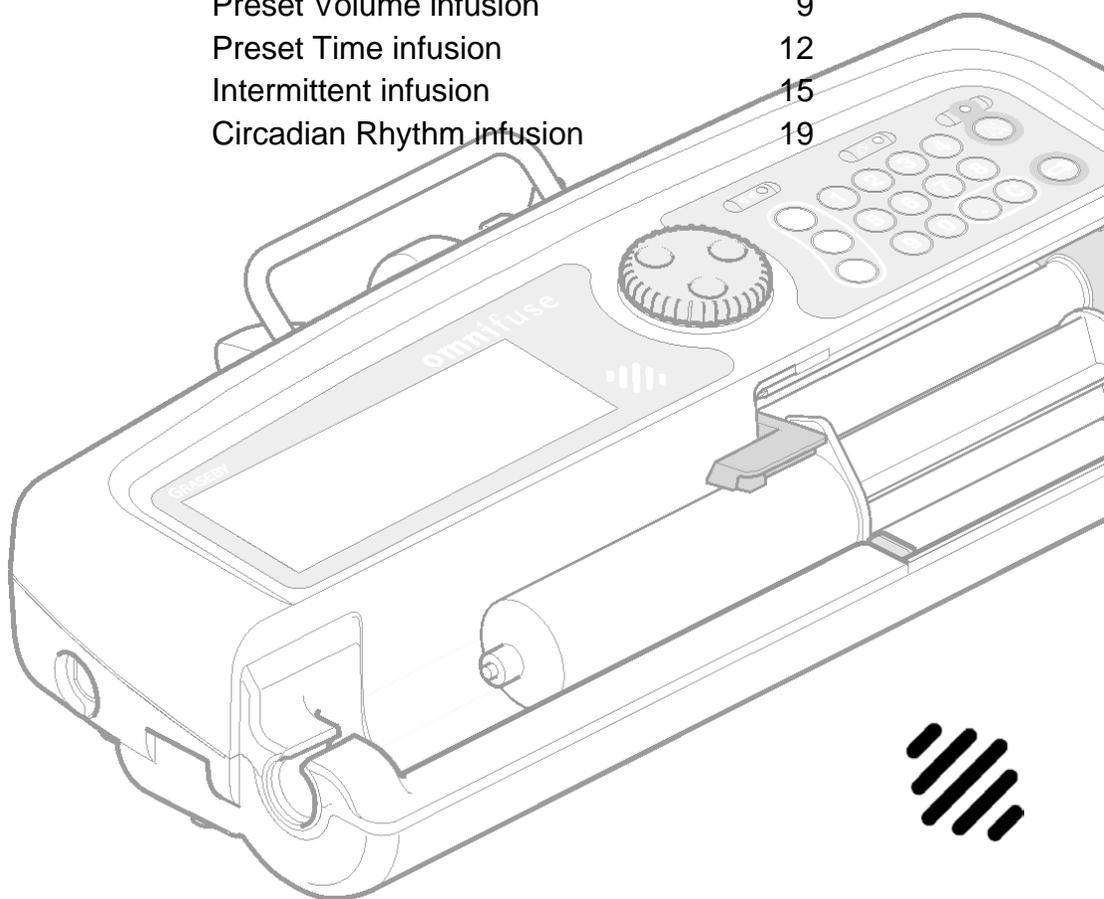

Mass units and special infusion modes

Continuous infusion in ml/h is always available on the Omnifuse pump.

Mass units

As well as programming in ml/h, the Omnifuse pump can also be programmed using mass units, for example mg/h or mg/kg/h.

To see a list of all the mass units that may be used on the Omnifuse, see the *Specification*, at the end of the manual.

Extra infusion modes

As well as the standard Continuous infusion, the Omnifuse pump can also be programmed for specialised infusions using the following infusion modes:

- Preset Volume
- Preset Time
- Intermittent
- Circadian Rhythm.

Setting up the mass units and extra modes

These specialised infusion modes and mass units will only be available for use on the Omnifuse if they have been set up by a technician.

The technician can do this using one of two methods:

- Configuration menu, for which a special access code and password are needed
- or
- Graseby Omnifuse Technician PC Software which can be used to set up the pump.

To find out about the Omnifuse Configuration menu and Technician software, see the *Omnifuse Technical User Manual*.

Programming using mass units

With the mass units option installed on the pump you can program the pump to reflect the way a patient's prescription has been written, taking into account how some drugs are dispensed.

- Drug mass per time period per kilogram of patient body-weight: for example mg/kg/h
or
- Drug mass per time period: for example µg/min.

Drug concentrations

Mass Unit infusions can be set up either with pre-mixed or with user-mixed drugs:

- Pre-Mix is when the drug is supplied already mixed with infusate solution in a specified concentration in mg/ml, µg/ml, or ng/ml
- User-Mix is when the drug and infusate solution are separate. You program the pump with the drug mass in mg, µg or ng, and the volume of the solution, and the pump then calculates the concentration for you.

The mass units feature can be used with Continuous, Intermittent and Circadian Rhythm modes.

Since the general principles of using mass units are the same for all infusion modes, this section illustrates their use by setting up a Continuous infusion.





Mass Units and the Totaliser

When you program an infusion using mass units, the totaliser will show you the total delivered ml and as a mass unit.

Total delivered since reset =	
2.2 ml (Bolus :	0.0 ml)
21.56mg (Bolus :	0.00mg)
Totals reset at 15:03 03/04/03	
THIS HOUR	HOLD

Setting up a Continuous infusion in mg/h

This example gives you an overview of how to set up a Continuous infusion using a pre-mixed syringe containing 10 mg/ml, infusing at 25 mg/h.

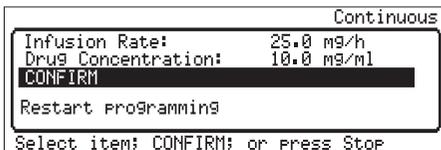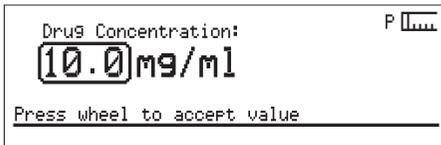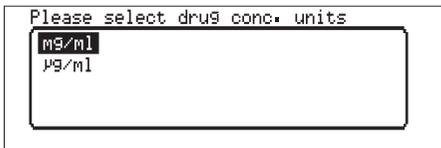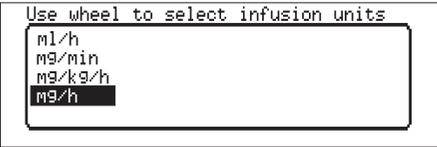

The instructions here start once you have loaded the syringe and confirmed the infusion mode (if necessary).

1. The screen shows a list of units available on this particular pump. Turn the Command wheel to select the correct units for this infusion, in this example, mg/h. Press the Command wheel.
2. On the next screen, enter the infusion rate according to the prescription, using the Command wheel or the numeric keypad. Press the Command wheel to confirm.
3. Select mg/ml as the drug concentration unit.
4. Enter the drug concentration and press the Command wheel to confirm.
5. When you have entered all the parameters, the Review screen is displayed.

Check the settings that you have entered, then **CONFIRM** to go to the Infusion Start screen.

Setting up a Continuous infusion in $\mu\text{g}/\text{kg}/\text{min}$

This example gives you an overview of how to set up a Continuous infusion at $2.5\mu\text{g}/\text{kg}/\text{min}$. The patient weighs 75 kg, and the concentration of drug in the syringe is 4 mg/ml.

The instructions here start once you have loaded the syringe and confirmed the infusion mode (if necessary).

Use wheel to select infusion units

ml/h
$\mu\text{g}/\text{kg}/\text{min}$
mg/h

Infusion Rate: P

2.5 $\mu\text{g}/\text{kg}/\text{min}$

Press wheel to accept value

Patient Weight: P

75.0 kg

Press wheel to accept value

Please select drug conc. units

mg/ml
$\mu\text{g}/\text{ml}$

Drug Concentration: P

4.0 mg/ml

Press wheel to accept value

Continuous




Infusion Rate:	2.5 $\mu\text{g}/\text{kg}/\text{min}$
Patient Weight:	75.0 kg
Drug Concentration:	4.00 mg/ml
CONFIRM	

Select item: CONFIRM: or Press Stop

1. The screen shows a list of units available on this particular pump. Turn the Command wheel to select the correct units for this infusion, in this case $\mu\text{g}/\text{kg}/\text{min}$. Press the Command wheel.
2. On the next screen, enter the infusion rate according to the prescription, using the Command wheel or the numeric keypad. Press the Command wheel to accept.
3. Enter the patient's weight and press the Command wheel to accept.
4. Select the drug concentration unit for the infusion - mg/ml and press the Command wheel to accept.
5. Enter the drug concentration and press the Command wheel to confirm.
6. When you have entered all the parameters, the Review screen is displayed. Check the settings that you have entered, then **CONFIRM** to go to the Infusion Start screen.


Setting up a Continuous infusion with user mix drug concentration

This example gives you an overview of how to set up a Continuous infusion when the drug and infusate solution are supplied separately.

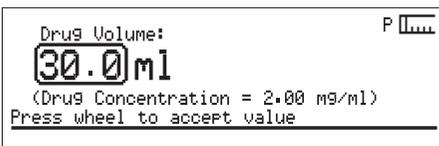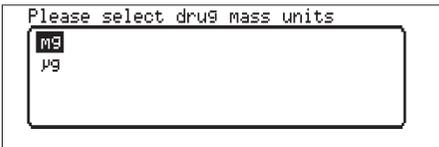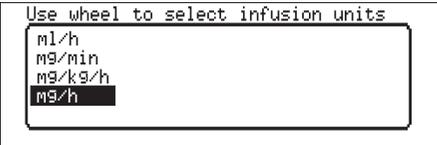

The instructions here start once you have loaded the syringe and confirmed the infusion mode (if necessary).

1. The screen shows a list of units available on this particular pump. Turn the Command wheel to select the correct units for this infusion, in this case mg/h. Press the Command wheel.
2. On the next screen, enter the infusion rate according to the prescription, using the Command wheel or the numeric keypad. Press the Command wheel to accept.
3. Select the drug mass units for the infusion - mg and press the Command wheel to accept.
4. Enter the drug mass and press the Command wheel to accept.
5. Enter the volume of solution in the syringe and press the Command wheel. As you enter the volume, the drug concentration is automatically calculated and displayed by the pump. Press the Command wheel to accept.




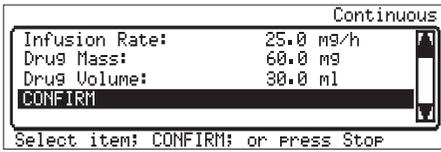


- When you have entered all the parameters, the Review screen is displayed. Check the settings that you have entered, then **CONFIRM** to go to the Infusion Start screen.

How to change the infusion mode

This explains how to change between the infusion modes available on the Omnifuse pump. The modes available on the pump shown here are: Continuous, Preset Volume and Preset Time.

If other modes have been set up on the pump, you may find that Intermittent or Circadian Rhythm infusion modes are also available on the menu. See *Setting up the mass units and extra modes*, page 4 - 1.

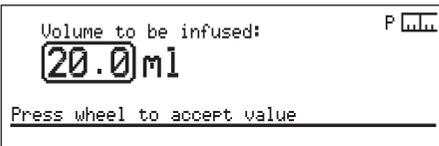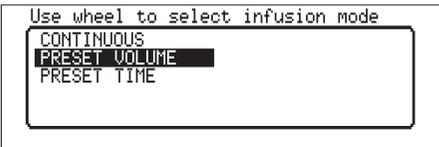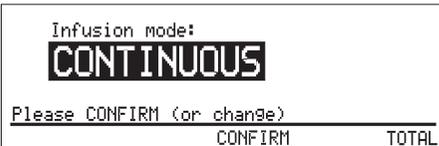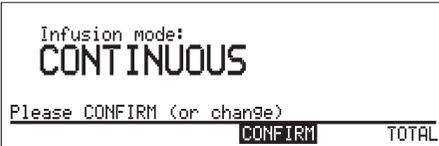

Where there is a choice of infusion modes, the first programming step is to confirm the mode, or change it.

1. To change it, turn the Command wheel to highlight the current infusion mode.
2. With the infusion mode highlighted, press the Command wheel.
3. A menu of infusion modes is displayed. Turn the Command wheel to highlight a new infusion mode.
4. Press the Command wheel to accept the new infusion mode.

The pump automatically displays the first programming screen for your chosen mode. This example shows the first screen for a Preset Volume infusion.


The Review screen

If you use an infusion mode with more than one parameter, you need to review and confirm your infusion settings.

The review screen will appear for Preset Volume, Preset Time, Intermittent and Circadian Rhythm infusions, and Continuous infusions if you are programming the pump using mass units.

Typical Review screen

A typical Review screen - in this case showing the parameters for a Preset Volume infusion - looks like this:

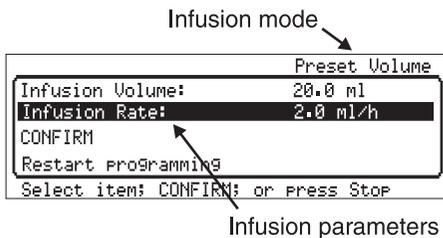


More complex Review screens

If there are more than five parameters, turn the Command wheel to review each one in turn. When you reach the end, you will be able to see the **CONFIRM** hotspot.

Viewing parameters during an infusion

You can check the settings when an infusion is running. Turn the Command wheel to highlight the **i** icon and press the wheel.

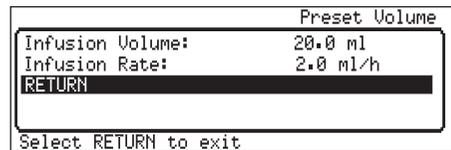

To exit from the information screen turn the Command wheel and select **RETURN**.

Using the Review screen

When the Review screen is displayed:


- You can correct any settings by highlighting the parameter and pressing the Command wheel to re-enter.
- If the settings are correct, turn the Command wheel to highlight **CONFIRM** and press the wheel. The Infusion Start screen will be displayed.

Restart programming

The Review screen also has a **Restart programming** hotspot, which you can use if you want to change the infusion mode.

Preset Volume infusion

This example shows how to program a Preset Volume infusion of 20 ml at 2 ml/h. On a pump which is used for Preset Volume infusions, you will see this screen once you have loaded the syringe:

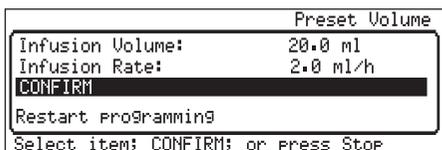

1. Press the Command wheel to **CONFIRM** the infusion mode. See page 4 - 7 if the pump displays a different infusion mode and you need to change to Preset Volume.
2. Program the Infusion Volume using the Command wheel or the numeric keypad.
Press the Command wheel to accept.
3. Enter the infusion rate and press the Command wheel to accept. The pump displays the Review screen so you may check the settings.
4. If either setting is wrong, highlight it and press the Command wheel to re-enter.
If the settings are correct, highlight **CONFIRM** and press the Command wheel.
5. The infusion is ready to start at this point, but you should first:
 - Check or change the Occlusion Alarm Level to ensure that it is appropriate for the infusion, see page 3 - 4.
 - Select **TOTAL**, to check or reset the cumulative total as shown on page 3 - 7.



You could also carry out these tasks if necessary:

- **REVIEW** the infusion settings.
- **PURGE** the line see page 3 - 10.
- Use the **SLEEP** feature see page 3 - 11.

6. To start the infusion press Start .

The running screen shows the infusion rate and the volume remaining to be infused.

Changing the rate while running - Preset Volume

Once the Preset Volume infusion is running, you can alter the infusion rate if necessary.



1. Highlight the current rate on the running screen and press the Command wheel.
2. Turn the Command wheel, or use the numeric keypad to enter the new rate.
Press the Command wheel to accept.
3. Press Start  to continue the infusion at the new rate.

If you do not press Start , the infusion continues at the previous rate, and the pump displays a warning message:

Rate change was not completed.

Checking the totals - Preset Volume

The message line shows the total volume infused since the last reset, including any amount delivered as a bolus.

To see the amount that was delivered as a bolus, highlight the **TOTAL** hotspot and press the Command wheel.

This shows the total since the last reset, and the bolus amount.

```
Total delivered since reset =
  20.0 ml (Bolus :   2.0 ml)

Totals reset at : 13:29 03/02/03
THIS HOUR THIS SHIFT  RESET  HOLD
```

If you want to see the hour-by-hour totals, highlight **THIS HOUR**.

For more details, see *The Totaliser*, page 3 - 6.

Giving a bolus - Preset Volume

The **BOLUS** hotspot appears once the infusion is running. Highlight the hotspot and press the Command wheel.

The bolus volume is deducted from the original volume to be infused.

For more details see *How to give a hands-on bolus*, page 3 - 15.

Ending a Preset Volume infusion



When the Preset Volume has been delivered an alarm sounds and a message is displayed that the infusion is complete.

Press Alarm Silence (A), then check the **TOTAL** and press and hold down Stop (C) to terminate the infusion.

```
INFUSION COMPLETE
The dose delivery is complete. The
infusion has been suspended.
```

Preset Time infusion

This example shows how to program a Preset Time infusion of 20 ml in 15 minutes. On a pump which is used for Preset Time infusions, you will see this screen once you have loaded the syringe:

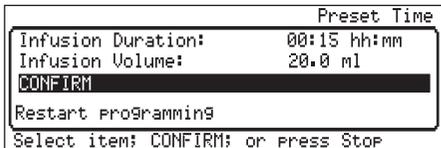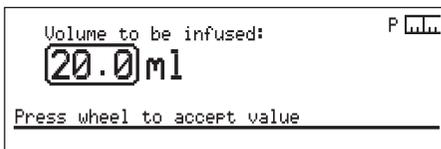

1. Press the Command wheel to **CONFIRM** the infusion mode.

See page 4 - 7 if the pump displays a different infusion mode and you need to change to Preset Time.

2. Enter the duration of the infusion using the Command wheel or the numeric keypad.

Press the Command wheel to accept.

3. Enter the infusion volume and press the Command wheel to accept. The pump displays the Review screen.

4. If either setting is wrong, highlight it and press the Command wheel to re-enter. If the settings are correct, highlight **CONFIRM** and press the Command wheel.

5. The infusion is ready to start at this point, but you should first:

- Check or change the Occlusion Alarm Level to ensure that it is appropriate for the infusion, see page 3 - 4.
- Select **TOTAL**, to check or reset the cumulative total, see page 3 - 7.



You could also carry out these tasks if necessary:

- **REVIEW** the infusion settings.
- **PURGE** the line see page 3 - 10.
- Use the **SLEEP** feature see page 3 - 11.

6. To start the infusion press Start .

The running screen shows the volume delivered so far.


Changing the rate while running - Preset Time

Once the infusion is running, you can alter the infusion rate by shortening or lengthening the duration of the infusion.




1. Turn the Command wheel to highlight the **CHANGE** hotspot and press the Command wheel.

The Duration screen appears, showing the original duration of the infusion.

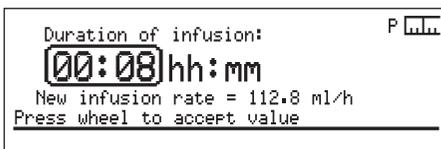

2. Turn the Command wheel, or use the numeric keypad to set the number of hours and minutes you want the infusion to run, from now. Press the Command wheel to accept the value.

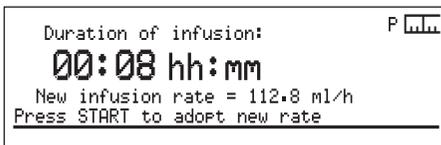

Check the new infusion rate based on the duration you have entered.



3. Press Start (⏪) to continue the infusion at the new rate.

If you do not press Start (⏪), the infusion continues at the previous rate, and the pump displays a warning message to remind you that the **Rate change was not completed**.

Giving a bolus - Preset Time

The **BOLUS** hotspot appears once the infusion is running. To give a bolus in addition to the original volume to be infused, highlight the hotspot and press the Command wheel.


The timer stops while the bolus is being delivered and restarts afterwards. The bolus volume is added to the volume delivered so far. For more details, see *Omnifuse bolus feature*, page 3 - 13.

Checking the totals - Preset Time

The message line shows the total volume infused since the last reset, including any amount delivered as a bolus.

To see the amount that was delivered as a bolus, highlight the **TOTAL** hotspot and press the Command wheel.

This shows the total since the last reset, and the bolus amount.

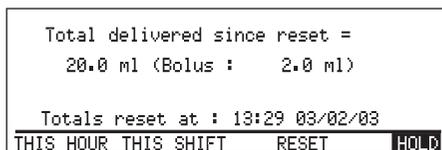

If you want to see the hour-by-hour totals, highlight **THIS HOUR**. For more details, see *The Totaliser*, page 3 - 6.

Ending a Preset Time infusion

When the infusion is finished, the pump displays this message and sounds an alarm.

- Press Alarm Silence (🔇), then press and hold Stop (⏹) to end the infusion.



Intermittent infusion


In Intermittent infusions, a dose is given at regular intervals over a defined period.

For example, the pump can be programmed to give a dose of 5 ml every four hours.

Optionally, you can program a start delay for an infusion.

Between doses the pump can carry on at a background rate, or wait without infusing until the next dose is due.

The diagram below shows an Intermittent infusion and explains the terminology involved.

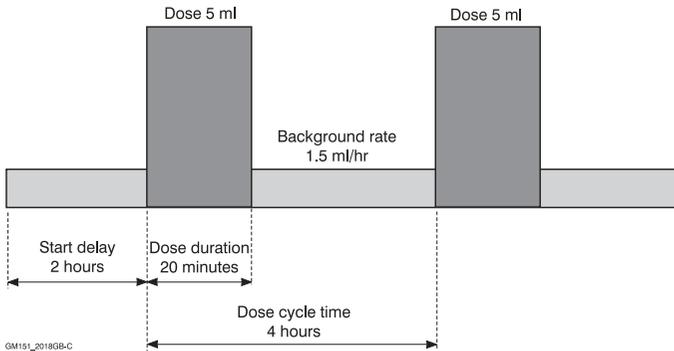

Elements of an intermittent infusion

Programming an Intermittent infusion

Infusion mode:
INTERMITTENT

Please CONFIRM (or change)

CONFIRM TOTAL

Dose: **20.0** ml

Press wheel to accept value

Dose Duration: **00:15** hh:mm

(Infusion Rate = 80.00 ml/h)

Press wheel to accept value

Dose Cycle Time: **01:00** hh:mm

Press wheel to accept value

Background Rate: **0.5** ml/h

Press wheel to accept value

Start Delay: **00:30** hh:mm

Press wheel to accept value

1. Switch on the pump and load the syringe, see *Basics, Chapter 2*.
2. **CONFIRM** the Intermittent infusion mode. If you need to change to Intermittent mode, see *How to change the infusion mode*, page 4 - 7.

3. Enter the dose and press the Command wheel to accept.

4. Enter the dose duration in hours and minutes. The Infusion Rate is calculated automatically. Press the Command wheel to accept.



5. Enter the dose cycle time in hours and minutes and press the Command wheel to accept.

6. Enter the background rate and press the Command wheel to accept. You can leave this set to 0.0 if no background infusion is required.

7. This screen appears if the pump has been set up to allow a Start Delay. Enter the start delay in hours and minutes and press the Command wheel to accept.

You can leave this set to 0.0 if no start delay is required for this infusion.

Starting an Intermittent infusion

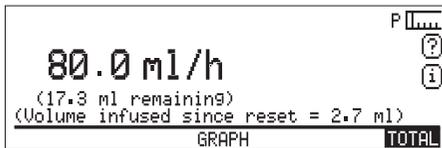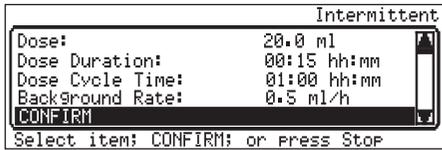

1. Check the settings for the infusion on the Review screen. Highlight **CONFIRM** and press the Command wheel to go to the Infusion Start screen.
2. On the Infusion Start screen, check the occlusion alarm level, see *How to set or change the occlusion alarm level*, page 3 - 4. You can also check or reset the **TOTAL** if required, see *The Totaliser*, page 3 - 6 and purge the line as described on page 3 - 10. To start the infusion press Start (⏪).

The pump begins to deliver the first dose. For pumps where a Start Delay has been set, see *Intermittent infusions with Start Delay*, page 4 - 18.

At the end of each dose, the pump displays a countdown to show the hh:mm until the next dose.

3. To terminate the infusion press and hold down Stop (⏹).

Infusion with Background Rate:

If a Background Rate has been programmed, the pump displays a screen like this after each dose to show the current infusion rate and time remaining until the next dose.



Note: If the pump has been set up so that a bolus is available, the **BOLUS** hotspot will be displayed.

Intermittent infusions with Start Delay



If you have programmed a Start Delay, when you press Start  the screen shows a countdown like this.

When the countdown reaches 00:00, the first dose is delivered.



If you have programmed a Start Delay **and** a Background Rate, the pump will begin the background infusion when you press Start .

The screen shows the background infusion rate and the countdown of hours and minutes shows how much time remains until the first dose will be delivered.

Background rates and KVO rate

If you have entered a background rate for your infusion, then depending on how the KVO rate has been configured by a technician, the pump will behave as follows when the syringe reaches the KVO point:

- If the KVO rate is higher than zero but *lower* than the background rate the pump goes into KVO rate until the syringe is empty
- If the KVO rate is *higher* than the background rate the pump will go into KVO mode, but will in fact infuse at background rate until the syringe is empty
- If the KVO rate has been set to zero the pump will continue at the current rate (dose or background), then stop and give a syringe-empty alarm.

Circadian Rhythm infusion

A Circadian Rhythm infusion is one where you program the pump to give a variable dose of medication each hour over a 24 hour period.

Unless you terminate the infusion it repeats itself the next day: that is to say, it finishes one cycle at 23:59 then starts the next one at 00:00.

This means that as long as the infusion is running and the syringe is replaced as required, the pump will go on repeating the 24-hour cycle indefinitely.

You can vary the dose between 0.1 ml and 99.9 ml every hour, or you can express each dose using mass units - see below. You can turn the infusion off for particular hours in the cycle by entering a rate of zero.

- The default dose is 1.0 ml/h.

When setting up a Circadian Rhythm infusion, the **VIEW LAST** and **COPY LAST** features allow you to copy the dose from one hour to the next. **VIEW LAST** also allows you to check the dose for the previous hour and correct it if necessary.


Circadian Rhythm and Mass Units

If required for the clinical application, the pump can be set up so you program a Circadian Rhythm infusion using mass units, for example mg/kg/h.

The screen shows you the dose that will be infused in the current hour of the day (for example 1 mg/kg) and also the rate calculated using the patient's weight (for example 8.5 ml/h for an 85 kg patient).




Setting up and running a Circadian Rhythm infusion

This example shows how to program the pump for a prescription which requires the drug to be given at 2.3 ml, increasing to 4.5 ml between 7.00 - 10.00.

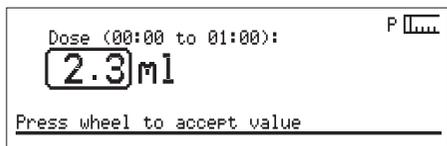

1. Switch on the pump and load the syringe, see *Basics, Chapter 2*.
2. **CONFIRM** the Circadian Rhythm infusion mode. If you need to change to Circadian Rhythm mode, see *How to change the infusion mode*, page 4 - 7.
3. The screen will display the first of the doses for the 24 hour period. Use the Command wheel to enter the dose of 2.3ml. Press the Command wheel to accept and go to the next hour, 01:00 to 02:00.
4. As the dose should be the same as the previous hour, you can use the **COPY LAST** hotspot rather than programming each dose. Highlight the hotspot and press the Command wheel to accept.
5. Repeat until 7:00-8:00, where you need to program the 4.5 ml dose. Program the 4.5 ml dose for 8:00 - 9:00, and 9:00 - 10:00.
6. Go through the remaining hours of the day programming the 2.3ml dose for each hour.

When you accept the final dose (for 23:00 to 24:00), pressing the Command wheel takes you to the Review screen.

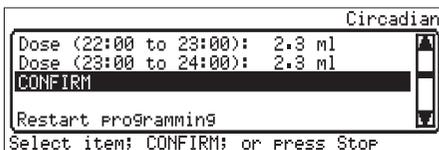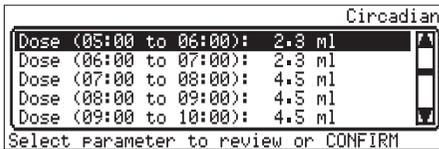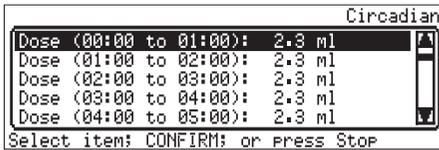

7. On the Review screen, you must view and confirm the doses for the whole 24 hour period before starting the infusion.

Initially, you can see the first five hours.

Turn the Command wheel to check the next five doses.

Continue turning the Command wheel and checking the remaining doses for the 24 hour period.

At the end of the list of doses, press the Command wheel to **CONFIRM**.

The screen displays the total to be infused over 24 hours.

If the prescription is correctly programmed, check or reset the **TOTAL**, then press Start (⏪) to start the infusion.

The screen shows the rate for the current hour.

8. To end the infusion press Stop (⏹).


If you do not end the infusion before 23.59 the pump will repeat the same cycle the next day.

Specialised infusions

Chapter 5

Drug Protocols

	<i>Page</i>
About Omnifuse Drug Protocols	1
Selecting a protocol	2

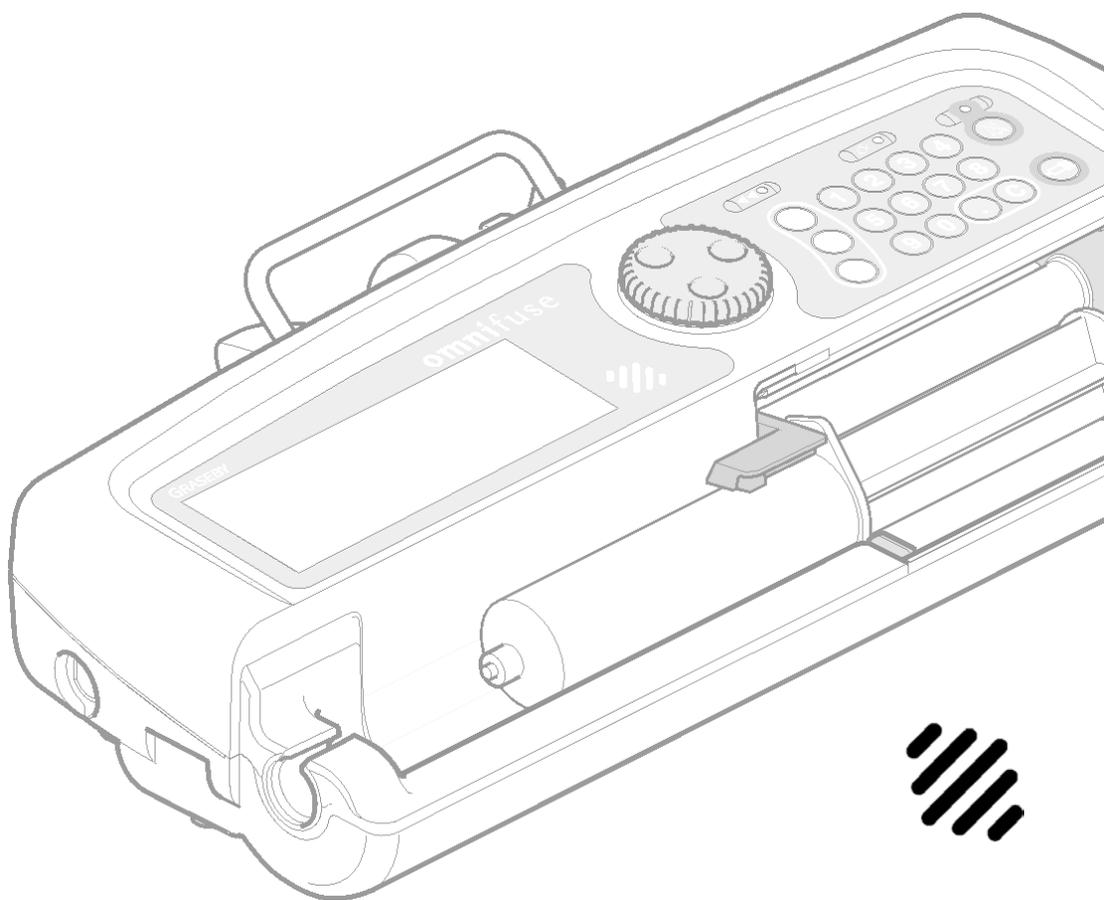

About Omnifuse Drug Protocols

The Omnifuse pump has a Drug Protocol feature which enables up to fifty predefined protocols to be stored in the pump. Each protocol specifies the prescription for administering a particular drug.

When selected on the pump, the settings for an infusion are presented ready for confirmation.

An example showing how you might select a protocol and run an infusion is shown on page 5 - 2.

Loading and storing a protocol on the Omnifuse

Before the pump is used in the clinical area, it is set up with the protocols that are required. The protocol must be set up and stored in the pump using the Omnifuse Drug Protocol Management System.

For details on how to load and store a protocol in the pump, see the *Omnifuse Technical User Manual*.

Programming using a protocol

This section gives an overview of how the pump operates when it has been set up to use protocols.

Select a protocol

The menu of protocols is displayed once you have loaded the syringe. The menu will contain only the protocols that have been set up for this pump. There may be just one, or up to fifty protocols.

Enter any required parameters

When you select a protocol, the pump displays either the Review screen, or the first parameter if the protocol requires input. Many of the parameters for an infusion will be set up as part of the protocol. The pump removes the majority of programming steps and only asks you to enter the minimum information it requires.

Review all the infusion parameters

When you have entered any parameters required by the protocol, the pump

displays all the infusion parameters.

This includes the parameters you have entered, and also the ones that are part of the protocol. You may be able to change some of these values, depending on how the protocol has been designed. Parameters in a protocol can be:

- Fixed, so you cannot make changes
- Limited, so you can make changes within a certain range specified in the protocol, or
- Changeable so you can alter them within the pump limits.

Start the infusion

The parameters must be confirmed before you can start an infusion.




When you press Start (⏪), the infusion begins and the pump displays the name of the protocol on the screen. The pump operates according to the confirmed protocol.

Selecting a protocol

The protocols available on each pump will have been set up for the clinical environment, so you can choose one that is suitable for your patient.

The screen shown here shows an *example* of how a previously named protocol may be used.

In this example, an infusion protocol called **Continuous 5 ml/h** has been set up in the pump for a clinical environment where it is used to administer a continuous infusion at 5 ml/h.

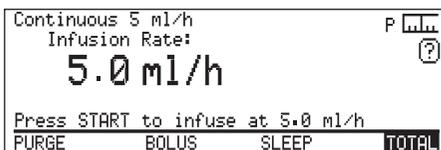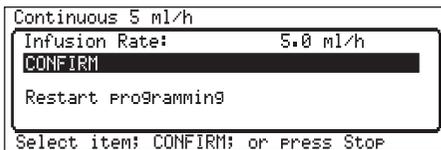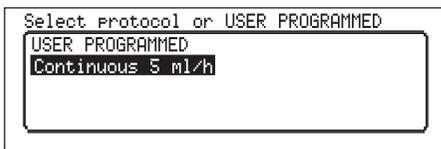

1. Switch on the pump and load the syringe, see *Basics, Chapter 2*.

As soon as the syringe is loaded, the pump shows a menu of one or more protocols.

2. Turn the Command wheel to highlight the name of the protocol you want to use. Press the Command wheel to select the protocol.

On the Review screen you must check that the infusion settings for the protocol match the prescription for the patient.

3. If correct, highlight **CONFIRM** and press the Command wheel.
4. Once you have confirmed the infusion parameters, the Infusion Start screen is displayed.

You should now carry out the usual pre-infusion checks described in *Programming the pump, Chapter 3*.

Note: Protocols are designed, named and loaded into the pump using the Omnifuse Drug Protocol Management System, see the *Omnifuse Technical User Manual*.

Programming the pump without a protocol

When Drug Protocols are in use, the pump presents a menu instead of the first programming step.

If a pump has been set up for use with Drug Protocols, but needs to be used for an infusion which is not covered by one of the protocols, it can be programmed with the **USER PROGRAMMED** option. This option appears on the Protocol menu.


USER PROGRAMMED can be used to set up a one-off infusion, by following the standard programming procedures described throughout this manual.

Note: The **USER PROGRAMMED** option can be turned off, so it does not appear on the Protocol menu. This is done using the Configuration menu, described in the *Omnifuse Technical User Manual*.

Chapter 6

In-line occlusion sensing

	<i>Page</i>
In-line Occlusion Sensing	1
Using the In-line pressure sensor	2
Monitoring In-line pressure	3

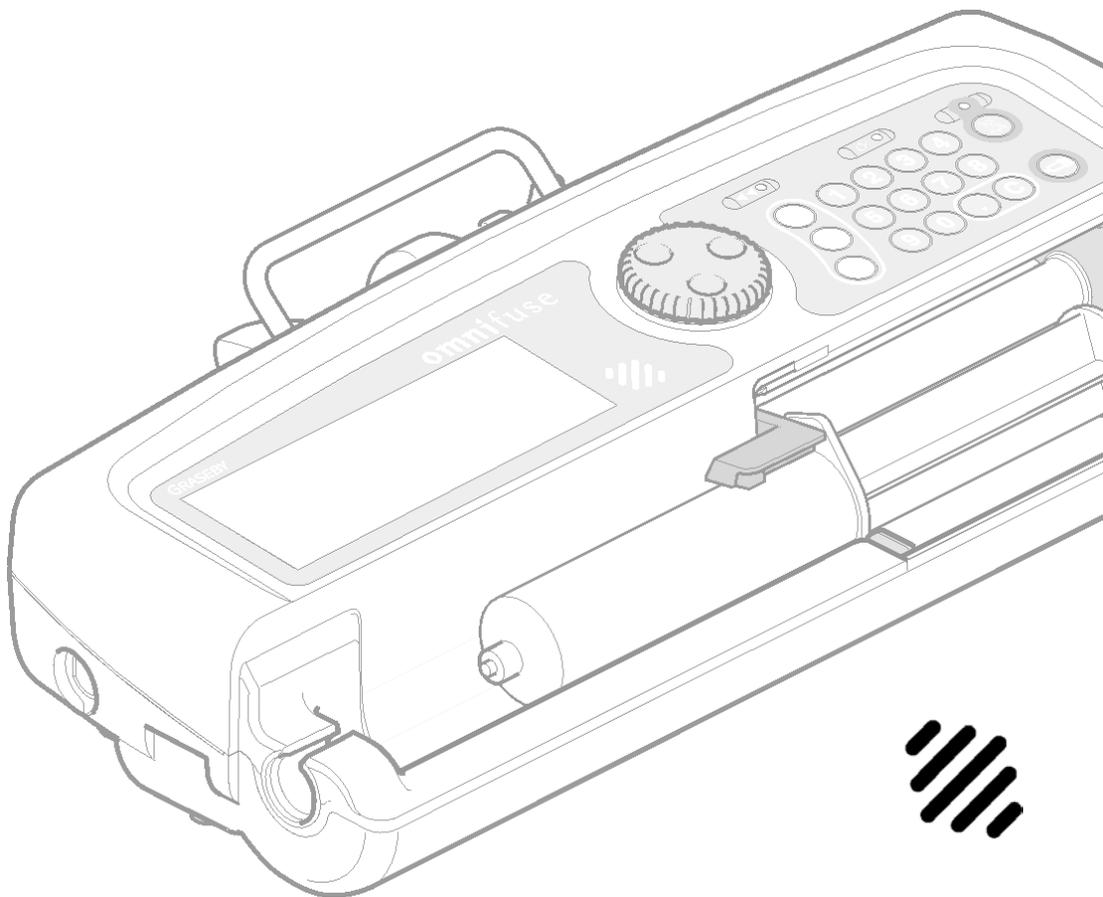

In-line Occlusion Sensing

The Omnifuse pump is available with an In-line occlusion sensor fitted.

On the standard Omnifuse pump, occlusion sensing works by measuring the back-pressure between the syringe plunger and the pusher block. This provides a good guide to the pressure level, as described in *Programming the pump, Chapter 3*.

Where faster warning of possible occlusions is required, for example a


neonatal unit or an ITU, the In-line Occlusion Sensor pump or In-line Sensor and Lockable cover pump should be used to detect pressure in the infusion line.

Special software in these pumps means that you can set the occlusion alarm to a specific value between 1 mmHg and 999 mmHg. If required, the upper limit can be restricted for a particular clinical area.

In-line pressure

To detect pressure in the line, you need to use a special extension set with a sensing disc which fits into the pump's pressure sensor.

The extension set is the Graseby Medical Flo- Safer™ extension set, part number 0130-0041:

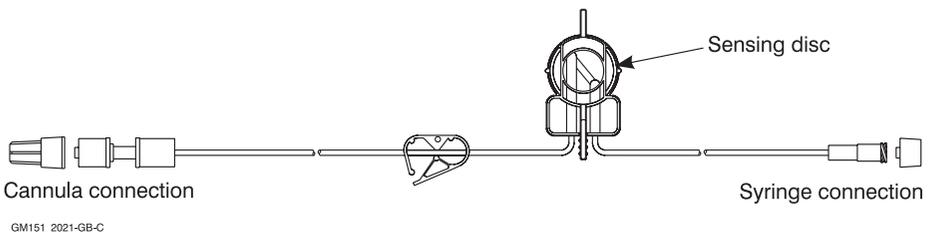


The set is suitable for use with a syringe of 5 ml capacity or larger. Always prime and connect the line according to the instructions for use on the Flo- Safer™ instruction leaflet.



WARNING: When using the In-line Occlusion Sensing option, use only a Graseby Medical or Graseby Medical - approved extension set (part number 0130-0041). Graseby Medical does not guarantee performance of the pump if unapproved extension sets are used. Failure to observe this warning may lead to compromised performance of the pump, resulting in patient injury or death.

Using the In-line pressure sensor

1. Ensure that there is no sensing disc in the In-line sensor, then switch on the pump.
2. Load the syringe and program the infusion following the usual procedure, see *Basics, Chapter 2*.
3. Insert the pressure sensing disc into the sensor. Press it firmly until you hear and feel the disc click into place.

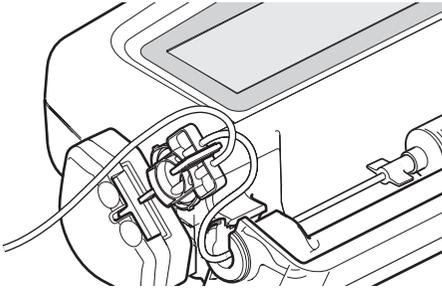

4. Set the occlusion alarm pressure: Turn the Command wheel to highlight the pressure icon then press the Command wheel.

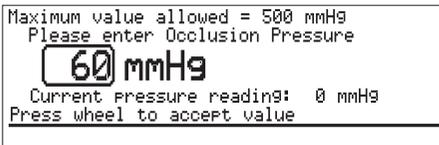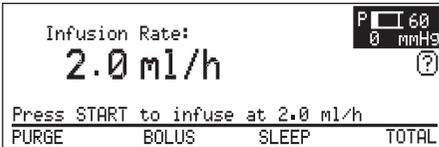


You will see a screen like this, showing the occlusion pressure alarm limit, the current alarm setting (boxed) and pressure reading.

5. Turn the Command wheel or use the keypad to enter the occlusion alarm pressure. Press the Command wheel to accept the alarm pressure.
6. Press Start (⏪) to begin the infusion.

Monitoring In-line pressure

The current in-line pressure and the alarm setting are displayed on the screen next to the In-line pressure icon:

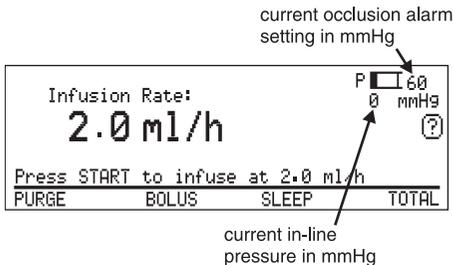

The In-line Occlusion Pressure icon shows how close the pressure is to the alarm level:

- If the pressure of fluid in the line is well below the alarm level, the icon will look like this: **P** 
- If the pressure is close to the alarm level and the Occlusion alarm is about to sound, the icon will look like this: **P** .

If the occlusion alarm pressure is reached, the pump will sound an alarm and suspend the infusion. The pusher will automatically move back from the syringe plunger to avoid accidental bolus delivery when the line is unclamped after clearing the occlusion.


Note: The Maximum value allowed for the Occlusion alarm is set on the Configuration menu by a technician, see the *Omnifuse Technical User Manual*.

Set the occlusion alarm pressure

To set the occlusion alarm pressure, turn the Command wheel to highlight the Occlusion Pressure display:



Press the Command wheel to display the Occlusion Pressure screen:

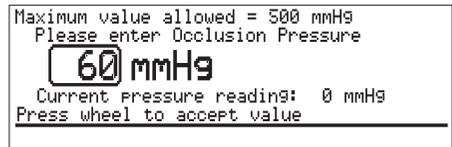

Turn the Command wheel or use the numeric keypad to alter the alarm setting, then press the Command wheel to exit the display.






Omnifuse with In-line Occlusion Sensing and Lockable Cover.

You can view, but not change the occlusion alarm setting if the infusion is running on an Omnifuse with In-line Occlusion Sensing and Lockable Cover .

You can change the setting if you first suspend the infusion, then unlock the cover.

How to view a Pressure Trend graph

If the Graphic feature is enabled on a pump which has In-line sensing in use, you will be able to view a Pressure Trend graph.

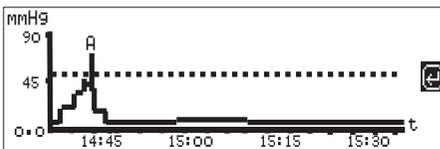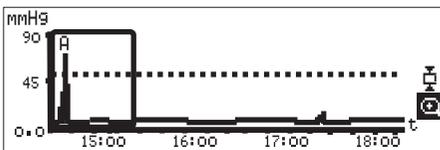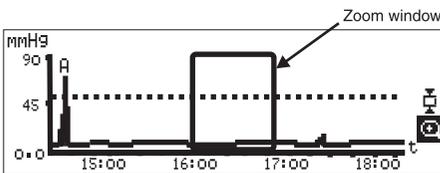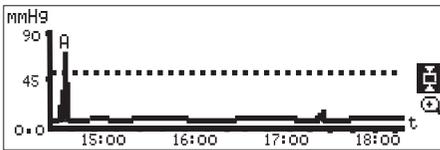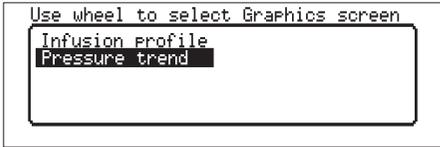

1. Select the **GRAPH** hotspot to see this menu. Highlight **Pressure trend** and press the Command wheel.
2. The Pressure Trend graph is displayed.
3. Select the Hold icon  to keep the graph on the screen for a minute. On the graph:
 - The dotted line shows the level at which the occlusion alarm has been set.
 - The unbroken line represents the pressure trend for the infusion.
4. Select the Zoom icon  to display the Zoom window.
5. Turn the Command wheel to move the Zoom window over the part of the infusion you wish to view.
6. Press the Command wheel to display the selected period in greater detail.
7. To leave the graph, press the Command wheel to select the Return icon .

In-line sensor warning messages

When using a pump which has an In-line occlusion sensor, you may see the warnings described here:

Message displayed:



Cause:

The disc is in the sensor when you switch on, which interferes with the pump’s self-tests.

Solution:

- Remove the pressure-sensing disc and press the Command wheel to continue.

If you repeatedly press the Command wheel without removing the disc, the pump will display a **DISC NOT REMOVED** message three times, then an **IN-LINE INVALID** message. The pump will use dry-side occlusion sensing.

Message displayed:



Cause:

In-line pressure sensing is enabled on the pump but you pressed Start  without inserting the pressure sensing disc.


Solution 1 - Continue without In-line sensing:

- Select **CONTINUE**. The infusion will begin, and for this infusion, the pump will operate as though the In-line sensor is not present.

Solution 2 - In-line sensing is required:

- Select **CANCEL**, or wait for the display to revert to the start screen
- Press the sensing disc firmly into the sensor until you hear and feel it click into place
- Press Start  to begin the infusion.

In-line occlusion sensing

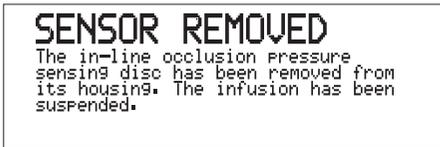In-line sensor warning messages, continued**Message displayed:****Cause:**

The sensing disc was removed while an infusion is running.

Solution:

- Press Alarm silence (A) to stop the alarm
- Replace the sensing disc
- Press Start (D) to restart the infusion.

Message displayed:**Cause:**

In-line pressure sensing is enabled on the pump, but you confirmed that the infusion should start without the In-line sensor. The pump displays this warning because it now detects that a sensing disc has been inserted.

Solution 1 - Continue without In-line sensing:

- Press Alarm silence (A) to stop the alarm
- Remove the sensing disc
- Press Start (D) to restart the infusion.

Solution 2 - In-line sensing is required:


You must end the current infusion and start again.

- Press and hold Stop (C) to end the infusion
- Insert the sensing disc
- Enter/Confirm the infusion parameters and press Start (D).

Chapter 7

Lockable cover

	<i>Page</i>
Introduction	1
About the security software	2
Removing and replacing the cover	5

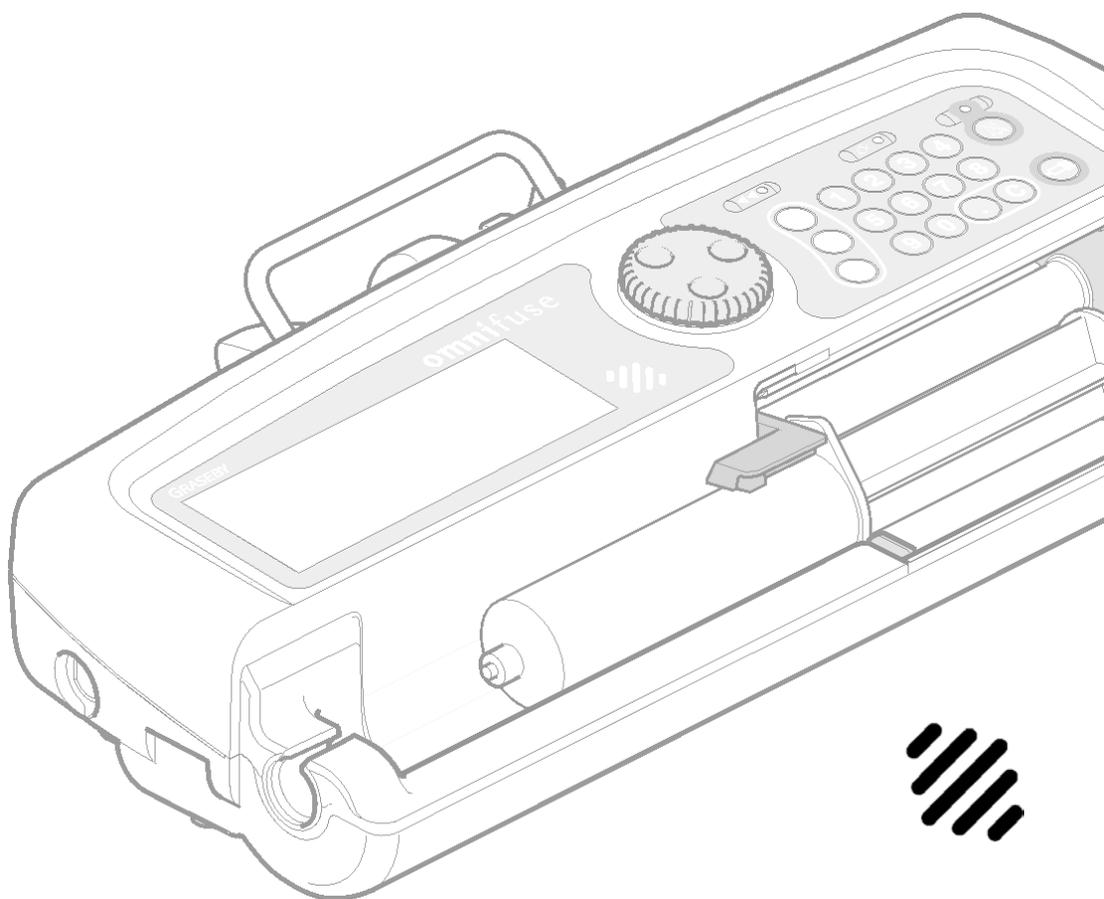

Introduction

This chapter describes how to use an Omnifuse pump with a Lockable cover. Both these Omnifuse pump variants are supplied with cover, lock and security software.

A standard Omnifuse pump may be upgraded with a Lockable cover, lock and security software.

Note: Diagrams in this chapter show the pump with Lockable cover only, although the instructions are the same for a pump with Lockable cover and In-line occlusion sensing. Contact your local Graseby Medical representative for further details.


Features of the Lockable cover

With a Lockable cover on the Omnifuse pump, the syringe can be secured in the pump during an infusion, as the cover can only be opened using a key.

Security software is provided with the Lockable cover. This can be used to prevent the infusion from starting until the cover is locked. The software also

locks the keypad once the infusion is started, so there can be no unauthorised change to the flow rate.

The security software can be switched off, even though the Lockable cover is fitted. In this case, you can use the cover to lock away the syringe but have a fully functioning keypad.

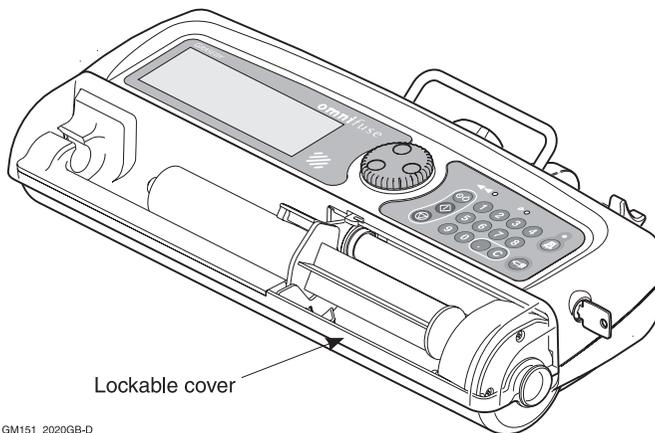

About the security software

Security software is supplied on all Omnifuse pumps with a Lockable cover.

Normally the security software and Lockable cover should be used together, as described on page 7 - 3.

However, if the Lockable cover on its own is secure enough for the clinical environment, the security software may be switched off by a technician, see the *Omnifuse Technical User Manual*.

Note: If you use a pump, with the security software switched off, the cover will secure the syringe, but the keypad remains active.

You can shut and lock the cover before or after starting the infusion, since the Start key (⏏) and the other keys on the keypad will operate normally.

The pump will not display any of the cover locked/unlocked warnings.

Is security software in use?


You can check whether or not the security software is in use on the pump as follows:

1. With the pump switched off and no syringe loaded, shut and lock the cover.
2. Switch on.
3. Check the message line. If the security software is in use, you will see the warning message: **Open cover to continue.**

WARNING: When the Omnifuse pump is fitted with a Lockable cover, the cover and associated security software should always be used in accordance with local protocol. The Lockable cover protects the syringe and its contents from tampering; the security software ensures that the pump settings are not interfered with. Failure to use the security software together with the Lockable cover could result in accidental or deliberate interference with pump settings, leading to patient injury or death.

Infusions with the cover *and* security software

The Lockable cover can be used with any of the Omnicuse infusion modes. When setting up an infusion on the Omnicuse with a Lockable cover and security software in use, follow these guidelines:

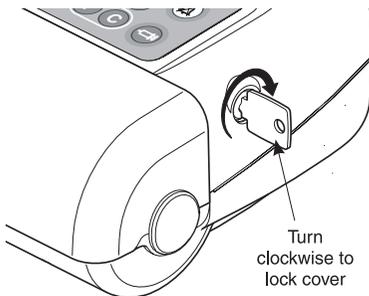

GM151_2022GB-D

1. When you load the syringe, ensure that the line leaves the pump through the channel, so it will not be trapped when you shut the cover.
2. With the cover open, follow the normal programming procedure for the infusion.
3. When the programming is complete, the message line reminds you to close the cover. However, before closing the cover:
 - Purge the line.
 - Check and change the occlusion alarm level.
 - Check and reset the totals.
 - Give a bolus if required.
4. Close the cover and turn the key to lock it shut. To lock, turn the key clockwise.

Always turn the key to ensure that the cover is locked shut. If you do not turn the key to lock the cover, you may be able to start the infusion, but the cover could fall open at any time. If this happens, the pump will stop infusing and sound the alarm.
5. Press Start (⏪) to begin the infusion.

COVER UNLOCKED

The infusion cannot be started with the cover open or unlocked.

6. If you have not locked the cover, the pump will remind you to do so.

Infusion running

Once the cover is locked and the infusion is running, you can:

- Suspend or stop the infusion
- Select the **TOTAL** hotspot to view the infusion totals
- Select the Information icon **i** to check the original infusion settings
- Silence the alarm.

Infusion suspended

You must suspend the infusion then unlock the cover if you need to:


- Change the flow rate
- Administer a bolus
- Change the occlusion alarm level.

Removing and replacing the cover

You should only remove the Lockable cover from the pump for cleaning. For details on cleaning the pump, see *Basics, Chapter 2, Care and maintenance, page 2 - 18.*

To remove the Lockable cover:

1. Unlock the cover and open it fully.
2. Load a syringe so that the pusher block is out of the way.
3. Pull the cover at the right-hand side until it flexes sufficiently to allow it to slide over the small boss on the end of the pump case.

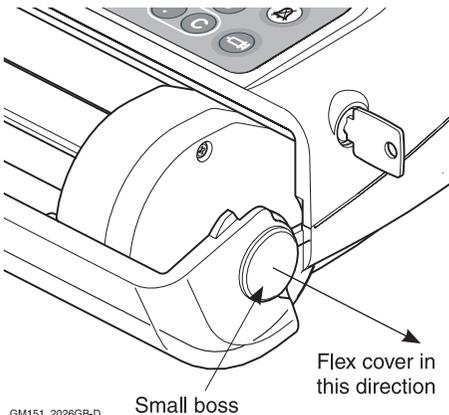

Initially you may find that you need to exert some force.

4. When the cover is free of the boss, pull it towards you then move it to the left to release it from the other end of the case.

To replace the Lockable cover after cleaning:


1. Locate the left-hand end of the cover in the hole in the left-hand end of the pump case, where the infusion line leaves the syringe trough.
2. Pull the right-hand end of the cover and flex it until it can slide over the small boss on the right hand end of the pump case.

Lockable cover

Chapter 8

Troubleshooting

	<i>Page</i>
Introduction	1
Alarms	4
Warnings	6
Repeated warnings	11

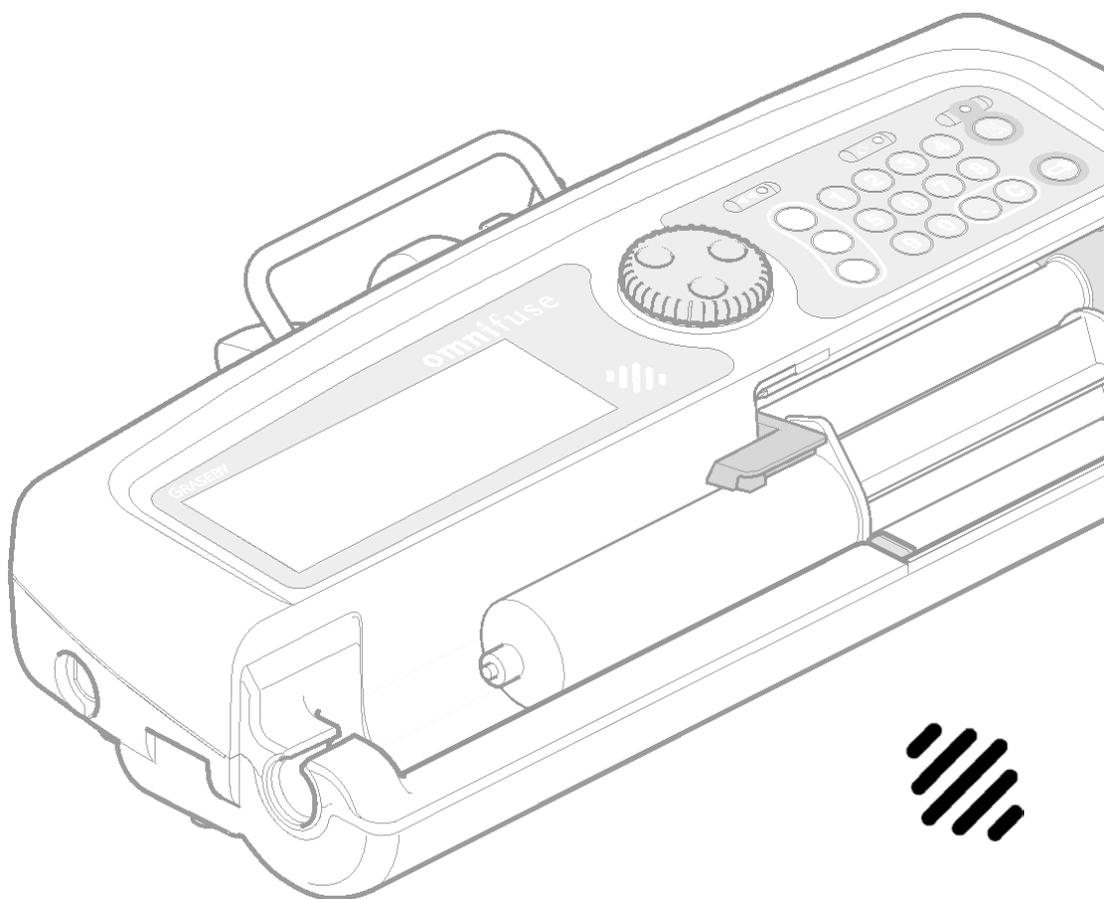

Introduction

This chapter offers troubleshooting advice if the pump sounds the alarm. It covers all the alarms and warnings that may be displayed on the Omnifuse screen to notify you of an operating or equipment problem.

The chapter also covers System Faults, which sound the alarm, but can only be resolved by a technician, see page 8 - 3.

Format of Omnifuse Alarms and Warnings

Omnifuse alarms and warnings are displayed on the screen in a standard format.

- The error message is displayed in large letters at the top of the screen
- Below this, on-screen information tells you how to clear the error.

This is an example of an alarm screen:

OBSTRUCTION

The syringe load/unload has been stopped as an obstruction has been detected. Please remove the obstruction and press LOAD to recommence the load/unload.

Lists of all alarms start on page 8 - 4 and warnings start on page 8 - 6.

How to handle Alarms

The pump sounds the alarm for conditions that require immediate attention. The LCD display flashes red and displays a message. If an infusion has been started, it is suspended.

Examples of conditions causing an alarm are:

- An operating problem, for example an occlusion has been detected, or the syringe is empty. This is the Occlusion alarm screen:

OCCLUSION

The infusion has been suspended. Please remove cause of occlusion and re-start.

- An external problem has stopped the infusion, for example someone has tampered with the syringe
- A preset infusion has ended. See a full list of Alarm messages on page 8 - 4.

If the pump sounds the alarm:

The screen always displays a message explaining what has caused the alarm.

1. Read the message.
2. Press Alarm Silence (Ⓜ) to silence the alarm and clear the message.
3. Check the patient and clear the problem which caused the alarm.
4. Restart the infusion.

How to handle Warnings

Warnings are displayed for conditions that require attention, but are not urgent problems. The alarm does not sound; instead the pump sounds a double beep. If the infusion has been started, it continues.


There are two types of warning, *repeated* and *one-off*, which are explained here. A list of warnings starts on page 8 - 6.

Repeated warnings


The pump displays the error message and on-screen information briefly, then shows the running screen, with the short message on the message line.

The error message and running screen alternate until you take action to resolve the warning condition.

For example, if you leave the pump running on battery power, and the batteries become depleted, you will see the warning message:

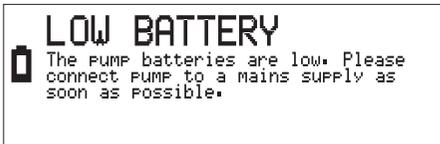

After three seconds this screen is replaced with the running screen where the message line displays a short message.

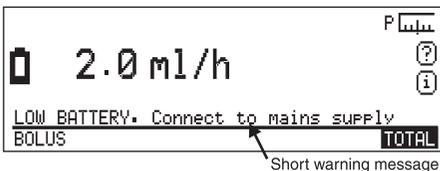


The message is repeated until you connect the pump to an AC mains power supply, or end the infusion and switch off.

One-off warnings

Warnings may be displayed while you are programming the infusion, or before the infusion is started, or when it is suspended.

These warnings are displayed on the screen for three seconds. You must follow the instructions on the screen before you can start the infusion.

For example, during the syringe load process, if the syringe ear is not correctly located the pump displays this message:

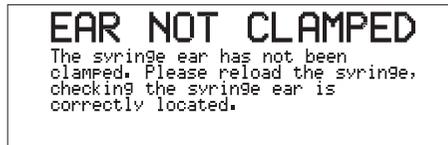

After three seconds this screen is replaced and the message line displays a short warning message:




If the pump displays a warning:

1. Read the message.
2. Carry out the action suggested on the screen.

How to handle System faults

The pump carries out a self-test when it is switched on and continually monitors all internal components while it is in use.

If it detects an error with one of the components the pump sounds a continuous alarm and displays a System Fault screen with a fault number and a message like the one shown here:

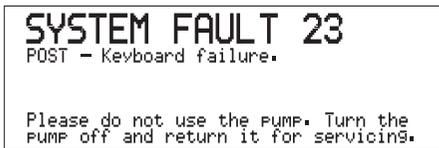

If the pump displays a system fault with a number:

1. The Alarm Silence key  has no effect. You must switch off the pump to silence the alarm. You do not need to note the fault number since this will be saved in the pump.
2. **Do not use the pump.** It must not be used again until it has been checked by a technician.

System fault at switch on

The system fault screen described opposite appears if a *component* fails at switch on or at any other time.

The following screen would be displayed only if a fault is discovered at switch on:



If you see this screen, you can press the Command wheel to display a help screen showing further information.

Do not use the pump. It must not be used again until it has been checked by a technician.

If the pump sounds a continuous alarm with no message:

When the pump is being used on battery, an independent battery monitor is active. This provides an additional check on the battery condition.

If the pump's alarm sounds but no message is displayed on the screen, the independent battery monitor has detected a problem.

The alarm cannot be silenced and the pump cannot be used. It must be taken to a qualified technician for repair.

Alarms

The pump sounds a continuous alarm and the screen flashes red when any of these messages are displayed.

You must press Alarm Silence  to clear the alarm. For more details, see *How to handle Alarms*, page 8 - 1.

Error message	On-screen information	Short message on screen
BATTERIES DEAD	Please connect pump to a mains supply immediately. The pump will automatically switch off.	URGENT: BATTERIES EXHAUSTED
CLAMP OPENED	The syringe clamp was opened during an infusion. The infusion has been suspended.	CLAMP OPENED
COMMS FAILURE	The remote connection to the pump has been broken.	COMMUNICATIONS FAILURE
CONTROL FINISHED	External controller terminated control whilst infusing.	CONTROL TERMINATED
COVER OPENED	The cover was opened during an infusion. The infusion has been suspended.	COVER OPENED
INFUSION COMPLETE	The dose delivery is complete. The infusion has been suspended.	Infusion complete, hold STOP key
LOAD NOT COMPLETE	The syringe loading operation did not complete within the time allowed. Open and close the barrel clamp then try loading the syringe again.	SYRINGE LOADING NOT COMPLETED IN TIME
OBSTRUCTION	The syringe load/unload has been stopped as an obstruction has been detected. Please remove the obstruction and press LOAD to continue the load/unload.	OBSTRUCTION
OCCLUSION	The infusion has been suspended. Please remove cause of occlusion and re-start.	OCCLUSION. See Help (?)

Error message	On-screen information	Short message on screen
SENSOR FITTED	The in-line pressure sensing disc has been fitted into its housing after the infusion has been started. The infusion has been suspended.	IN-LINE OCCLUSION SENSOR FITTED
SENSOR REMOVED	The in-line occlusion pressure sensing disc has been removed from its housing. The infusion has been suspended.	IN-LINE OCCLUSION SENSOR REMOVED
SYRINGE EMPTY	The infusion has been suspended. Please load a new syringe or end this infusion.	SYRINGE EMPTY. See Help (?)
	The preset dose has not yet been delivered. Please change the syringe and press START to deliver the remaining dose.	DOSE INCOMPLETE, LOAD NEW SYRINGE
	The infusion is incomplete. Please change the syringe and press START to deliver the remaining dose.	INFUSION INCOMPLETE, LOAD NEW SYRINGE
SYRINGE TAMPERING	The syringe has been disturbed. The infusion has been suspended.	SYRINGE TAMPERING DETECTED
	The syringe has been disturbed. Open and close the barrel clamp then try loading the syringe again.	

Warnings

There are two lists of warnings in this section: One-off warnings, and Repeated warnings.

For more information these types of warning and details on what to do if the pump displays a warning, see *How to handle Warnings*, page 8 - 2.

One-off warnings

Error message	On-screen information	Action
ALARM STATE	The previous alarm has not been cleared. Clear alarm condition, press ALARM SILENCE , or press STOP as appropriate.	When this warning message disappears from the screen, read the short message line. It will tell you the alarm state that needs to be cleared.
BRAND LOCKED	The pump can only be used with this brand of syringe.	If you want to load a syringe of a brand other than the one displayed, ask a technician to change the pump configuration.
CONC. INVALID	The entered or derived drug concentration is too high/low. The screen indicates the max/min concentration allowed.	Check the settings you have programmed against the patient's prescription.
CONFIRM / REVIEW	The infusion programming sequence has not been completed. Ensure that all infusion parameters have been confirmed.	Select REVIEW to check/change the infusion parameters and CONFIRM them. You can then press Start  .
COVER LOCKED	This option cannot be selected when the cover is locked.	Press Stop  to suspend the infusion before unlocking the cover.
COVER UNLOCKED	The infusion cannot be started with the cover open or unlocked.	Shut the cover and lock to secure it before pressing Start  .
DELAY INVALID	The start delay is too high. The maximum allowed is __	Enter a shorter delay (in hours:mins) for the start of the intermittent infusion.

Error message	On-screen information	Action
DISC NOT REMOVED	Please remove the disc from in-line sensor. In-line calibration cannot be carried out without removing the disc.	The pump was switched on with an in-line sensing disc in the sensor. Remove it before pressing the command wheel to continue.
DOSE INVALID	The infusion dose, bolus dose or dose mass is too high/low. The maximum/minimum dose allowed is _	Check the settings you have programmed against the patient's prescription.
DOSE START	A pre-programmed dose has started. The bolus request has been aborted.	If a bolus is required, wait until the dose has been delivered and try again.
DRUG MASS INVALID	The drug mass is too high/low. The Maximum/minimum drug mass allowed is _	Check the settings you have programmed against the patient's prescription.
DRUG VOL. INVALID	The drug volume is too high/low. The Maximum/minimum drug volume allowed is _	Check the settings you have programmed against the patient's prescription.
DURATION INVALID	The infusion, dose or bolus duration is too short/long. The minimum/maximum duration allowed is _	Check the settings you have programmed against the patient's prescription.
EAR NOT CLAMPED	The syringe ear has not been clamped. Please reload the syringe, checking the syringe ear is correctly located.	Load the syringe again and place the syringe ear, or flange, in the syringe ear slot.
FIXED PROTOCOL	This parameter (or the whole protocol) is fixed and cannot be altered.	You cannot change the selected parameter in this protocol. Accept this protocol/parameter as it is, select a different protocol, or choose USER PROGRAMMED .
HISTORY DOWNLOAD	The infusion cannot be started whilst a history download is in progress.	Wait for the download to end, or press stop if the pump should not be downloading history.

Error message	On-screen information	Action
IN-LINE INVALID	In-line sensor calibration cannot be carried out and therefore is invalid. Dryside sensor measurement will be used instead.	If you need to use in-line occlusion sensing, remove the disc from the sensor, switch off the pump then switch on again.
INVALID DATE/TIME	The start/stop date or time you entered is incorrect. Please correct.	Re-enter the date and time for the History display.
INVALID RATE	The bolus rate is invalid. Please adjust the bolus time or dose.	Check the settings you have entered for the preset bolus. Re-enter the duration/size of dose, as the settings you have entered resulted in a rate that is too fast.
INVALID SYRINGE	Please replace the syringe.	You must use a syringe of the same brand if you need to replace a syringe during an infusion.
INVALID TIME	The bolus duration is invalid. Please adjust the bolus dose or rate.	Check the settings you have entered for the preset bolus. Re-enter the infusion rate/size of dose, as the settings you have entered resulted in a duration that is too short or long.
KVO RATE REDUCED	The configured KVO rate is greater than the infusion rate or the background rate. The KVO rate has been reduced to an appropriate value.	For your information only, you do not need to take any action as the automatic adjustment means the infusion will start.
NO SYRINGE	None of the standard syringe brands are enabled. At least one should be enabled. Select and change at least one brand to ENABLED.	Only appears when a technician configures the pump for use, see <i>Omnifuse Technical User Manual</i> .
NO SYRINGE LOADED	BOLUS option was selected without a syringe loaded. Please load a syringe before starting.	Check that the syringe is loaded before you select BOLUS .

Error message	On-screen information	Action
NO RECORDS FOUND	No history records exist between the dates specified. Please select alternative dates, or exit history.	Check the range of dates you have entered.
PASSWORD CHANGE	The Configuration password has been changed. The Technician password has been changed.	Only appear when the technician is configuring the pump, see the <i>Omnifuse Technical User Manual</i> .
PASSWORD ERROR	The entered password is the same as one of the other passwords. All passwords must be different. Please enter another password.	
PASSWORD MISMATCH	Entered passwords do not match. The password has not been changed.	
PROTOCOL ERROR	An error occurred whilst trying to load the selected protocol. Please select another protocol, or select USER PROGRAMMED if that is available.	You should consider having the pump checked by a technician since all protocols available should load successfully.
RATE REDUCED	The default bolus rate has been reduced to the maximum allowed by the loaded syringe and configured maximum bolus rate.	The maximum bolus rate is restricted by the size of the loaded syringe. If the maximum bolus rate is not suitable for the prescription, you must use a larger syringe.

Error message	On-screen information	Action
RATE INVALID	The background rate entered is too high/low for the loaded syringe. The maximum/minimum rate allowed is _____	Check the setting you have made for the background rate, infusion rate or bolus rate against the patient's prescription.
	The derived or entered infusion rate is too high/low. The screen indicates the max/min rate allowed.	
	The bolus rate is too high/low. The maximum/minimum allowed is _	
	The new infusion rate is lower than the KVO rate. The minimum infusion rate is _	
SYRINGE LOADED	Syringe brands cannot be enabled or disabled whilst a syringe is loaded. Unload syringe first.	Only appears when a technician configures the pump for use, see <i>Omnifuse Technical User Manual</i> .
TIME INVALID	The dose duration must be less than the dose cycle time.	Check the setting you have made for the intermittent dose duration/cycle time against the patient's prescription.
	The dose cycle time is too high/low. The maximum/minimum dose cycle time allowed is _	
	The entered time is invalid. Only times in the range 00:00 to 24:00 may be entered.	Enter the correct time using the 24-hour clock. Remember to enter 4:00 as 04:00.
VOLUME INVALID	The infusion volume is too high/low. The maximum/minimum volume allowed is _____	Check the setting for the preset infusion volume against the patient's prescription.
WEIGHT INVALID	The patient weight is too high/low. The Maximum/minimum patient weight allowed is _	Check the patient's weight against the weight you have entered.


Repeated warnings

Error message	On-screen information	Short message on screen
_ min TO END	The syringe is nearly empty. There is now less than _ minute(s) to the end of the infusion.	LESS THAN __ min(s) TO END
_ ml TO END	The syringe is nearly empty. There is now less than _ ml to the end of the infusion.	LESS THAN __ ml TO END
AC MAINS FAILURE	The mains power supply has been disconnected. The pump is now running on battery.	AC MAINS FAILURE
BACKUP BATTERY	The backup battery should be replaced. The current date/time may be lost.	BACKUP BATTERY LOW
BATTERY FAULTY	The pump cannot be used on battery.	BATTERY FAULT. Do not use on battery
COMMS FAILURE	The remote connection to the pump has been broken.	COMMUNICATIONS FAILURE
KVO	Infusion ended. Pump now running in KVO mode at _	NOW RUNNING AT KVO RATE __
LOW BATTERY	The pump batteries are low. Please connect pump to a mains supply as soon as possible.	LOW BATTERY. Connect to mains supply
NOT INFUSING	The pump is not infusing. Please set up and start the infusion or switch the pump off.	NOT INFUSING
PUMP TOO COLD	The internal temperature of the pump is lower than expected. Please move pump to a warmer environment.	PUMP TOO COLD
PUMP TOO HOT	The internal temperature of the pump is higher than expected. Please move pump to a cooler environment.	PUMP TOO HOT

Error message	On-screen information	Short message on screen
RATE NOT CHANGED	The START key was not pressed. The infusion rate has not been changed. Enter new rate again then press START.	RATE CHANGE NOT COMPLETED
SOUNDER FAILURE	No sound can be detected. Remove pump from service as soon as possible.	SOUNDER FAILURE
VOLUME LOW	There is insufficient fluid remaining in the syringe to deliver another complete bolus dose.	VOLUME INSUFFICIENT FOR NEXT BOLUS

Specification/Standards

	<i>Page</i>
Specification	1
Standards	8

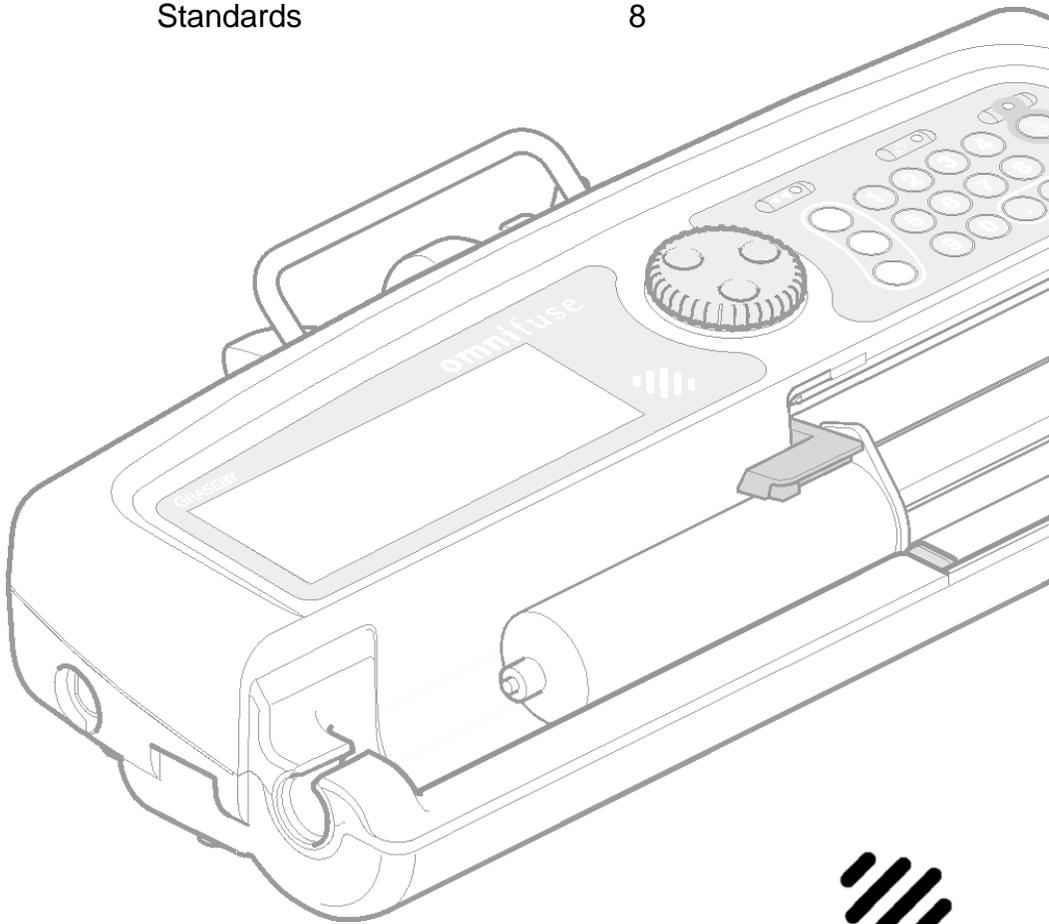

Specification

Weight	3.5 kg (approx)
Dimensions	384 mm x 170 mm x 92 mm (not including pole clamp)
Orientation	Horizontal, either mounted on a pole or flat on a stable horizontal surface
Display	LCD super twist with viewable area of approximately 105 mm x 32 mm. Green backlight when connected to AC mains supply and for up to 3 minutes (configurable) following a key-press when operating on battery. Backlight flashes red during an alarm unless configured Off
Data retention time	More than 12 months
Operating temperature	5° to 40°C
Storage temperature	-20° to 55°C
Relative humidity	20% to 90% non-condensing (operating)
Atmospheric pressure	700 - 1060 millibars (operating)
Alarm volume	Louder than 65 dBA at 1 metre at maximum volume
Software options	Extra Infusion modes and Mass Units Drug Protocols Graphics Remote Control

Pump variants

Description	Part numbers
Omnifuse Syringe Pump:	
UK and Ireland Model	0159-0001
Australian Model	0159-0740
Canadian (English) Model	0159-0711
Omnifuse with In-line Occlusion Sensing:	
UK and Ireland Model	0152-0001
Australian Model	0152-0740
Canadian (English) Model	0152-0711
Omnifuse with Lockable Cover:	
UK and Ireland Model	0157-0001
Australian Model	0157-0740
Canadian (English) Model	0157-0711
Omnifuse with In-line Occlusion Sensing and Lockable Cover:	
UK and Ireland Model	0158-0001
Australian Model	0158-0740
Canadian (English) Model	0158-0711

Power supply

AC power supply	100 - 240 V at 50/60 Hz. 50 W
Battery type	Set of three Cyclon™ sealed lead-acid batteries
Battery operating time	10 hours at 5 ml/h
Battery charge time	10 hours
Backup battery	Single 3 V lithium battery

Infusion flow rates

Range	0.1 to 800 ml/h dependent on syringe size
Accuracy	±2% measured over the 2nd hour of an infusion at 1 ml/h and at 5 ml/h with a Braun Omnifix 50 ml syringe and 150 cm extension set
Bolus accuracy	±5% with a Braun Omnifix 50 ml syringe and 150 cm extension set measured over 25 boluses of 1 ml each
KVO rate	Between 0.05 ml/h and 2 ml/h
Purge rate	50, 100, 200, 400, 800 ml/h (upper limit dependent on syringe size)
Bolus rate	0.1 to 800 ml/h (upper limit dependent on syringe size) in increments of: <ul style="list-style-type: none"> • 0.1 ml/h up to 100 ml/h • 1 ml/h above 100 ml/h

Maximum infusion flow rates

Syringe Size (ml)	Flow Rate (ml/h)
2	50 ml/h
3	50 ml/h
5	100 ml/h
10	200 ml/h
20	400 ml/h
25	400 ml/h
30	600 ml/h
50/60	800 ml/h

Programming ranges

Continuous infusion rate programming ranges by infusion unit

Unit	Range	Increment
ml/h	0.1 to 800	0.1 up to 100, 1 above 100
ml/min	0.01 to 13.0	0.01 up to 10, 0.1 above 10
mg/kg/h	0.1 to 99.9	0.1
µg/kg/h	1 to 999	1
ng/kg/h	1 to 999	1
mg/h	0.1 to 500	0.1 up to 100, 1 above 100
µg/h	0.1 to 999	0.1 up to 100, 1 above 100
ng/h	0.1 to 999	0.1 up to 100, 1 above 100
mg/kg/min	0.01 to 99.9	0.01 up to 10, 0.1 above 10
µg/kg/min	0.1 to 999	0.1 up to 100, 1 above 100
ng/kg/min	1 to 999	1
mg/min	0.01 to 50.0	0.01 up to 10, 0.1 above 10
µg/min	0.01 to 99.9	0.01 up to 10, 0.1 above 10
ng/min	1.00 to 99.9	0.01 up to 10, 0.1 above 10

Infusion pressure

Maximum infusion pressure 1250 mmHg

infusion pressure

Note: This value is approximate and is the pressure at the front face of the syringe plunger. It also assumes an ideal syringe with no stiction and a low infusion rate.

Occlusion sensing - pressure levels

Alarm levels	Approximate pressure in mmHg
Level 1	180
Level 2	300
Level 3	500
Level 4	750
Level 5	1250

Note: Approximate pressure at the front face of the syringe plunger assuming an ideal syringe with no stiction and a low infusion rate.

Occlusion sensing - time to occlusion

Alarm level and infusion rate	Number of minutes to occlusion	Bolus on occlusion release
Level 1 @ 1 ml/h	15	Less than 0.1 ml
Level 1 @ 5 ml/h	2.5	
Level 5 @ 1 ml/h	More than 100	Less than 1.0 ml
Level 5 @ 5 ml/h	21	

Note: Values are approximate and are determined using the method described in EN 60601-2-24 clause 51.6b using a Braun Omnifix 50 ml syringe and a 150 cm line (part number 0128-0122).

Note: Values are for a 50 ml syringe. Values are reduced for smaller syringes.

Accessories

Flo-Safer™ extension sets	Length	Part number
Syringe extension sets	150 cm	0128-0122
	200 cm	0128-0198
Syringe extension sets with anti-syphon valve	150 cm	0128-0253
	200 cm	0128-0254
Syringe extension set with pressure sensing disc	150 cm	0130-0041
Polyethylene-lined syringe extension sets with anti-syphon valve	150 cm	0128-0257
	200 cm	0128-0258
Low-priming volume syringe extension set with anti-syphon valve	100 cm	0128-0259
Epidural syringe extension sets with anti-syphon valve (Yellow)	150 cm	0128-0261
	200 cm	0128-0262
Epidural syringe extension sets (Yellow)	150 cm	0128-0263
	200 cm	0128-0264

Software	Part number	
Graseby Omnifuse Technician PC software	0151-0266	
Omnifuse Drug Protocol Management System	0153-0084	
Omnistack	Part number	
Omnifuse pump stacking system	0156-0001	
Wheelbase	0156-0096	
Pole assembly for wheelbase	Long	0156-0097
	Short	0156-0098

Supported syringe brands and sizes

Brand/Size	2	3	5	10	20	30	50/60
BD Plastipak							
BD Precise							
Braun Euroject							
Braun Omnifix							
Braun Perfusor							
Codan							
Faulding Pharmaject*							
Fresenius Injectomat							
IMS Pumpjet*							
JMS							
Monoject							
Nipro							
Terumo							
TOP							

* Pre-filled syringe

Note: The syringes shown in the **Supported syringe brands and sizes** table are supported by Omnifuse with their critical dimensions but may not achieve the stated accuracy due to syringe variability (with the exception of the Braun Omnifix 50ml).




Symbols

Symbol	Meaning
	The applied part is Type CF
	Identifies the Potential Equalisation Terminal located on the body of the pole clamp
	The pump should be operated from an AC power source

Standards

	Electrical safety	Classified as Internally Powered Equipment: Class 1, Type CF insulation on all inputs
Fluid ingress protection		IPX4 Splash-proof
	CE marking	The CE mark demonstrates that the pump conforms to the requirements of European Council Directive 93/42/EEC concerning medical devices. The number 0473 identifies the Notified Body under which the Quality Systems operated within Graseby Medical Ltd. are assessed
Design standards		EN60601-1, EN60601-1-2, EN60601-1-4, EN60601-2-24
	Disposal	<p>When the time comes to dispose of the pump, its batteries or any of its accessories, do so in the best way to minimise any negative impact on the environment. You may be able to use recycling or disposal schemes. To find out about these, contact your local waste disposal service.</p> <p>Separate any other parts of the equipment where arrangements can be made for their recovery; either by recycling or energy recovery.</p> <p>The only pump components which are potentially harmful enough to require separate disposal according to manufacturer's instructions or local regulations are:</p> <ul style="list-style-type: none"> • Main batteries (lead acid) • Back-up battery on main PCB (lithium) • LCD display (contains harmful chemicals and may explode if incinerated) <p>Note: Existing national or local regulations concerning waste disposal must take precedence over the above advice</p>
Patents		Applied for

Startup curves

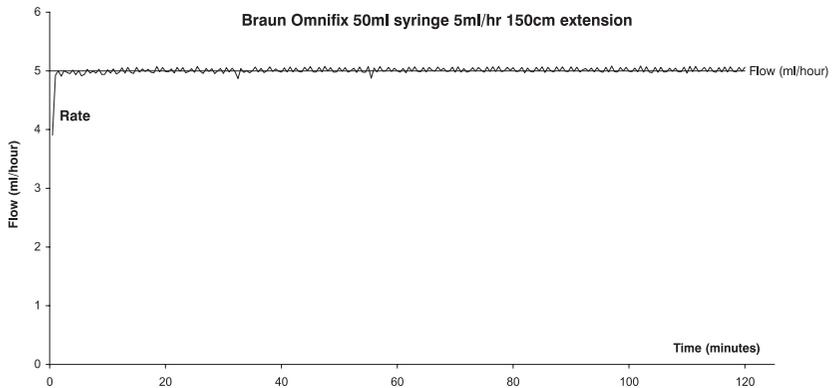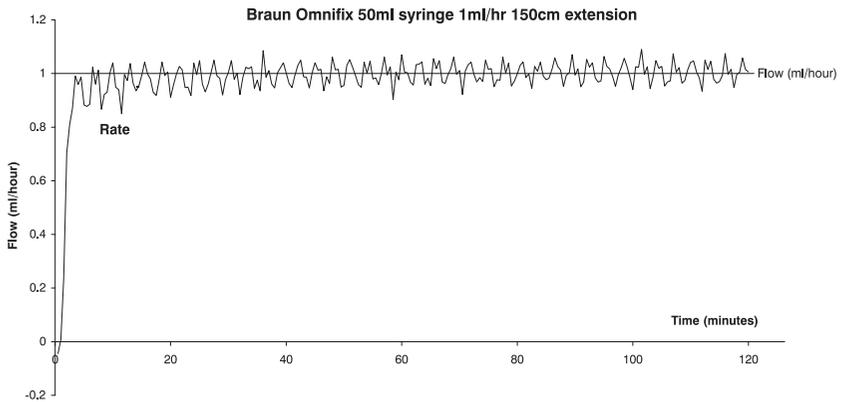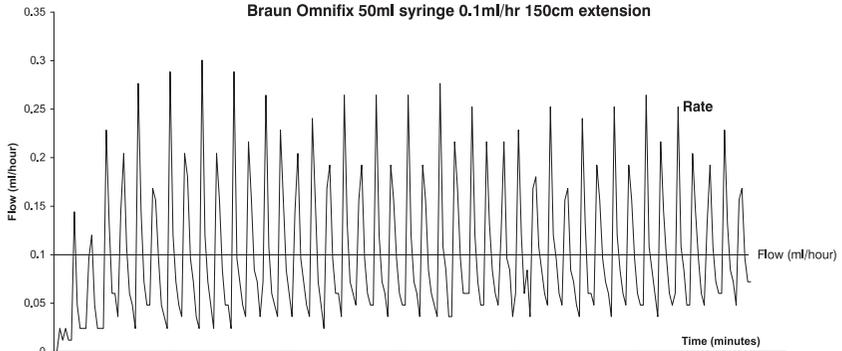


GM151_2102GB-A


Specification

Trumpet curves


Braun Omnifix 50ml syringe 0.1ml/hr 150cm extension

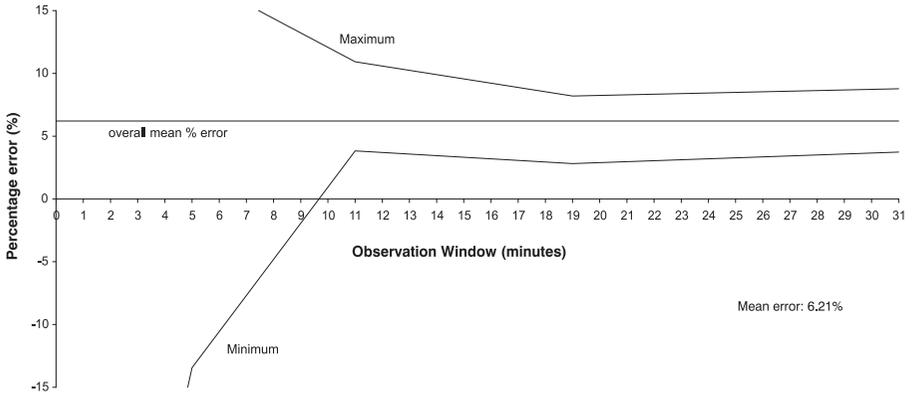

Braun Omnifix 50ml syringe 1ml/hr 150cm extension

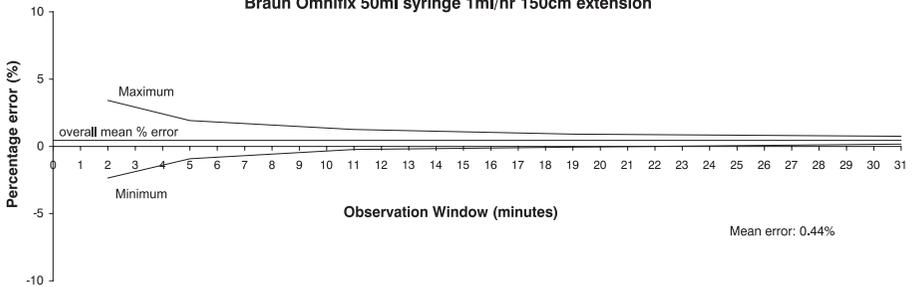

Braun Omnifix 50ml syringe 5ml/hr 150cm extension

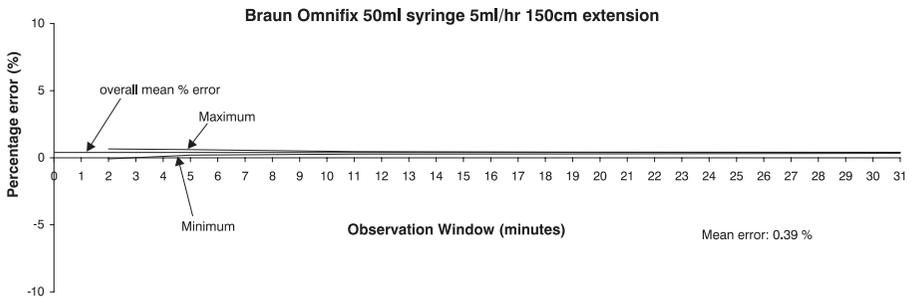

GM151_210IGB-A

Specification

Index

A

- AC power supply S-2
- Accessories S-6
- Administration sets 2-13
- Alarm
 - infusion complete 4-11, 4-14
 - LED 2-3, 2-10
 - occlusion level setting 3-4
 - occlusion pressure 3-3
 - Silence key 2-3
 - sounds 2-11
 - volume 2-11
- Asleep
 - pump 3-11

B

- Battery
 - charging time S-2
 - disposal S-8
 - icon 2-9
 - level 2-7
 - recharging 2-6
 - type S-2
- Battery life 2-7
- Beeps
 - adjusting the level 2-11
- Before 2-4
- Bolus

- continuous infusion 4-9, 4-12
 - graph 3-17
 - hands-on 3-15
 - in preset infusions 4-11
 - occlusion pressure during 3-3
 - occlusion release 6-1, S-5
 - overview 3-13
 - pre-set 3-16
 - rate range S-3
 - totals 3-6, 3-14
- Bolus features 1-3

C

- Carrying the pump 2-18
- CE marking S-8

Change

- brand of syringe 2-17

Change syringe

- during an infusion 2-13

Circadian Rhythm infusion

- main steps 4-20

Clamp pump to pole 2-18

Cleaning the pump 2-20

Command wheel and keypad 2-3

Concentration, drug 4-2

- pre-mixed 4-3
- user-mixed 4-5

Configuration 1-2

Continuous infusion

- programming (mass units) 4-3

Cover

- removing 7-5

D

Delayed start to infusion 4-18

Design standards S-8

Dimensions of pump S-1

Disposal of batteries S-8

Dose

- circadian 4-20
- intermittent 4-16
- mass units 4-2

Drug protocols 5-1

E

Electrical safety S-8

Embolism, avoiding 3-10

External communications 1-3

F

Faulty pump message 2-4

Flo-Safer S-6

Fluid ingress S-8

G

Graph

- icons 2-9, 3-18
- infusion profile 3-19
- menu 3-17

pressure trend 6-4
text info 3-19
zooming 6-4

H

handle 2-18
Handle, carrying 2-18
Help icon 2-9, 3-18
Hotspots 2-8
Hourly total 3-8

I

Icons 2-9
Indicators 2-10
Info icon 2-9, 3-18, 3-19
Infusion

- background rate 4-17
- change infusion mode 4-7
- change rate 3-1
- change rate when running 3-2
- circadian rhythm mode 4-19
- extra modes 4-1
- indicator 2-10
- intermittent mode 4-15
- optional modes 4-1
- viewing totals 3-7

In-line pressure

- monitoring 6-3
- overview 6-1
- sensing disc 6-2
- setting alarm pressure 6-3

IPX4 S-8

IV Pole 2-18

K

KVO rate 4-18

L

LCD S-1
LEDs 2-3, 2-10
Load syringe 2-15
Lockable cover

- optional 7-1

M

Mains

- indicator LED 2-10
- lead retainer 2-6
- supply 2-6
- supply failure 2-6

Mains supply S-2
Maintenance 2-20
Mass units

- drug concentration 4-2
- examples 4-3, 4-4
- programming ranges S-4

Maximum infusion rate S-3
Message line 2-8
Monitoring

- battery 2-6
- occlusion pressure 6-3

N

Near empty

- infusion indicator 2-10

O

Occlusion

- alarm 2-11
- alarm level 3-4
- dry-side 3-3
- icon 2-9
- in-line 6-1
- in-line alarm level 6-2

Omnistack 2-18
Options

- Extra Infusions 4-1, 5-1, 6-1
- graphs 3-17
- in-line occlusion sensing 6-1
- lockable cover 7-1
- protocols 5-1

P

Parts of the pump 2-1
PC software overview 1-5
Pole clamp

- standard 2-18

Power

- battery 2-6
- power-up tests 2-4

- summary 1-3
- supply S-2
- switch on pump 2-4

Preset

- bolus 3-16

Previous hour total 3-9

Protocol

- saving 5-2
- selecting 5-2

Pump configuration 1-2

Pump ID 2-4

Purge

- how to purge line 3-10
- rate range S-3
- volume 3-10

R

Rate

- change 3-1
- change when running 3-2

Recharging batteries 2-6

Remove syringe 2-16

Review

- complex screens 4-8
- program settings 4-8
- settings on WAKE ME 3-12

S

Safety features 1-1

Saving a protocol 5-2

Service due message 2-4

Shift total 3-8

Size

- of pump S-1
- of syringes S-6

Sleep mode 3-11

Software

- overview 1-5
- security 7-2

Sound

- adjusting the level 2-11
- alarms and warnings 2-11

Standards S-8

Startup curves S-9

Switching off 2-5

Switching on 2-4

Syringe

- brand, changing 2-17
- brands supported S-6
- change shown on graph 3-18
- how to load 2-15, 2-16
- Load key 2-3
- removing when switched off 2-16
- replacing during an infusion 2-13
- securing with cover 7-1
- size/rate/increments S-3

T

Titrate infusion rate 3-2

Total

- bolus amount 3-14
- cumulative 3-6
- hour by hour 3-8
- hourly 3-8
- previous hour 3-9
- shift 3-8
- view or reset 3-6

Trumpet S-10

Trumpet curves S-10

U

User-mixed drug concentration 4-2

V

Volume

- alarm 2-11
- alarm adjusting 2-11
- preset infusion 4-9

W

Wake me 3-12


Warnings

- one-off 8-6
- repeated 8-2

Weight

- patient's 4-4
- pump S-1

For further information, please contact your local distributor or
Graseby Medical direct on +44 (0)1923 246434



GRASEBY

Graseby Medical Ltd.

Colonial Way, Watford, Herts, UK, WD24 4LG

Telephone: +44 (0)1923 246434, Facsimile: +44 (0)1923 231595

<http://www.smiths-medical.com>

CE
0473

Part No. 0151-0839-C

January 2004

© 2004 Graseby Medical Ltd

smiths

Smiths Medical - A part of Smiths Group plc