SERVICE MANUAL

SERIES 160

Order no.

350-8108	Index b	Service manual cpl.	Service-Anleitung kpl.
350-7901	Index b	Content only	nur Inhalt

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Subject to technical modification

Copyright

These service instructions are protected by copyright. All rights reserved. No part of the service instructions may be reproduced, processed, duplicated or published in any form by photocopying, microfilming, reprinting or other process, in particular electronic means, without the written agreement of Precisa Instruments AG.

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Section E: Adjustment - - - - - E1

1. Open up the balance	E2
2. Removal of the base plate	E3
3. Replacing the beam.....	E4
3.1 Remove the beam.....	E4
3.2 Installing the beam	E4
3.3 Adjustment and check of over- and underload.....	E4
4. Inspection / Adjustment.....	E5
5. Check and adjustment of over-, and underload protection.....	E6
5.1 Checking overload	E6
5.2 Checking under-load	E6
6. Check corner load	E7
6.1 Adjustment of corner load	E7
6.2 Example of corner load adjustment.....	E7
6.3 Linearity of the balance	E9
6.4 Final check of the adjustment	E9
7. Programming the Flash-memory with the balance specific Software nnnn-nnn.PEE.....	E10
7.1 Installing the EEPROM writer	E10
7.2 Programming the Flash-memory.....	E10

8. Data sheet - adjustment and final control	E11
9. Programming the Flash-memory with the balance program software snn-nnnn.BIN	E12
9.1 Installing the Precisa „DownloadVnnn.exe“ program	E12
9.2 Programing the Flash-Memory.....	E12
10. Error messages	E13
10.1 Operating error (Error disappears with correct use or application)	E13
10.2 Fatal error (Balance stays still, error repair in service mode).....	E13
10.3 Hardware error (Balance stays still)	E13
11. Discussion with the manufacturer	E14

Section A: Accessories

Kapitel A: Zubehör

1. Accessories Serie 160 / Zubehör Serie 160


Note: The accessories are identical as for the model listed in the „Service Manual Accessories“

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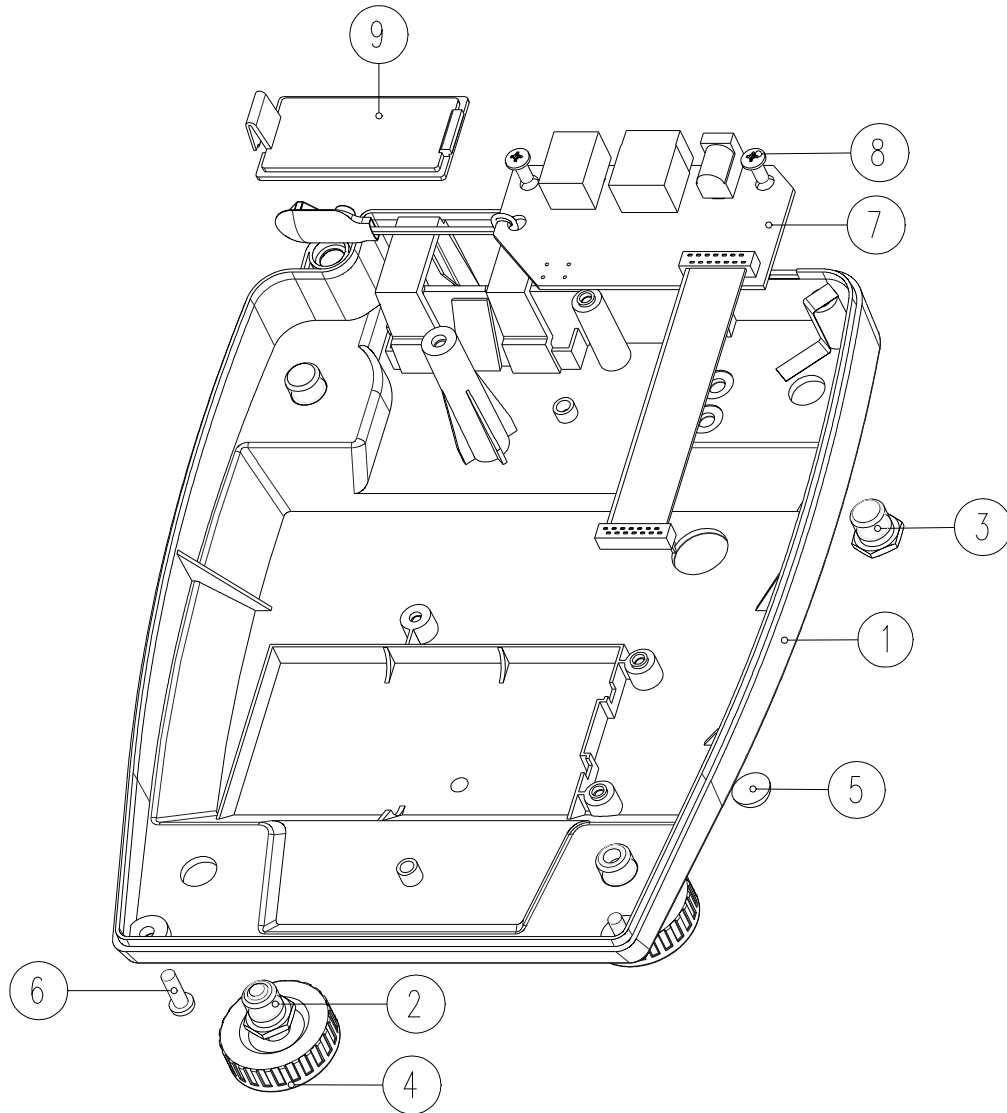
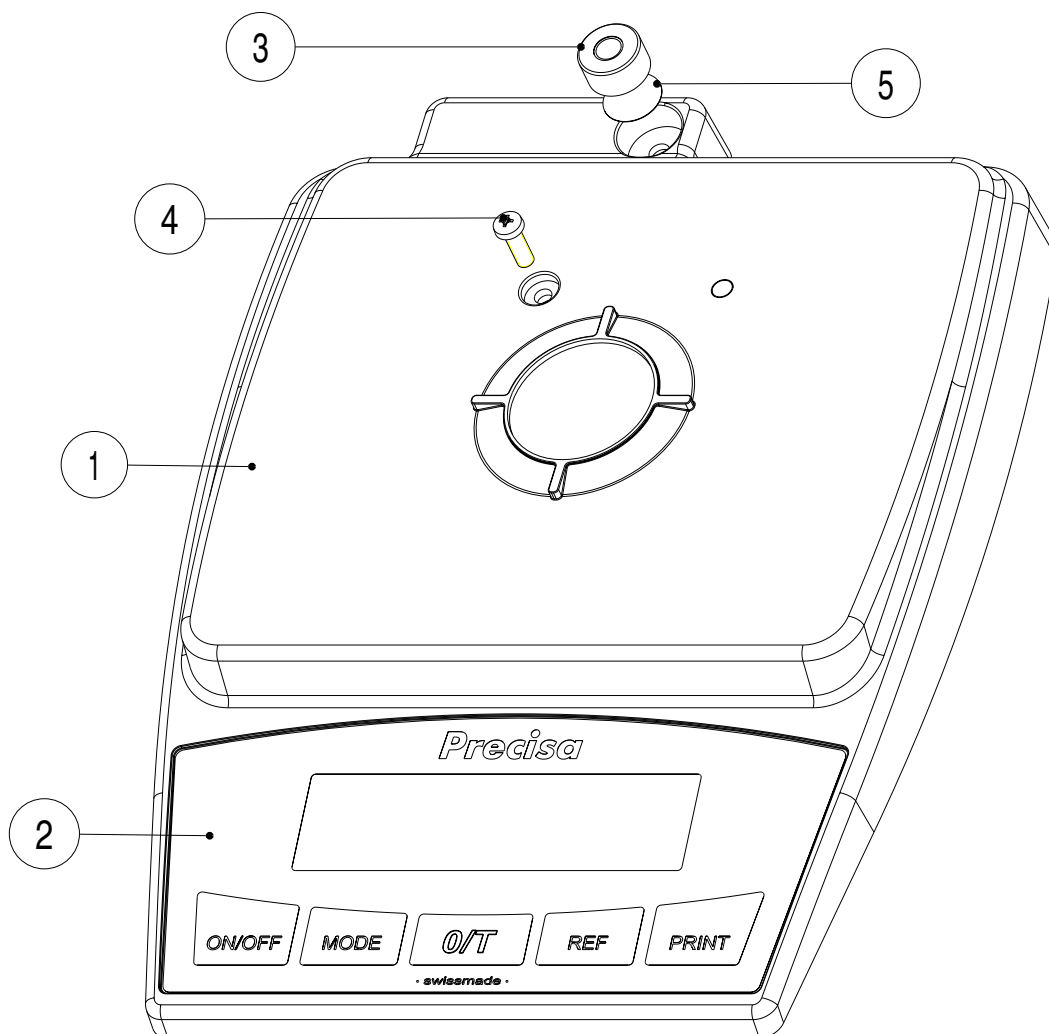
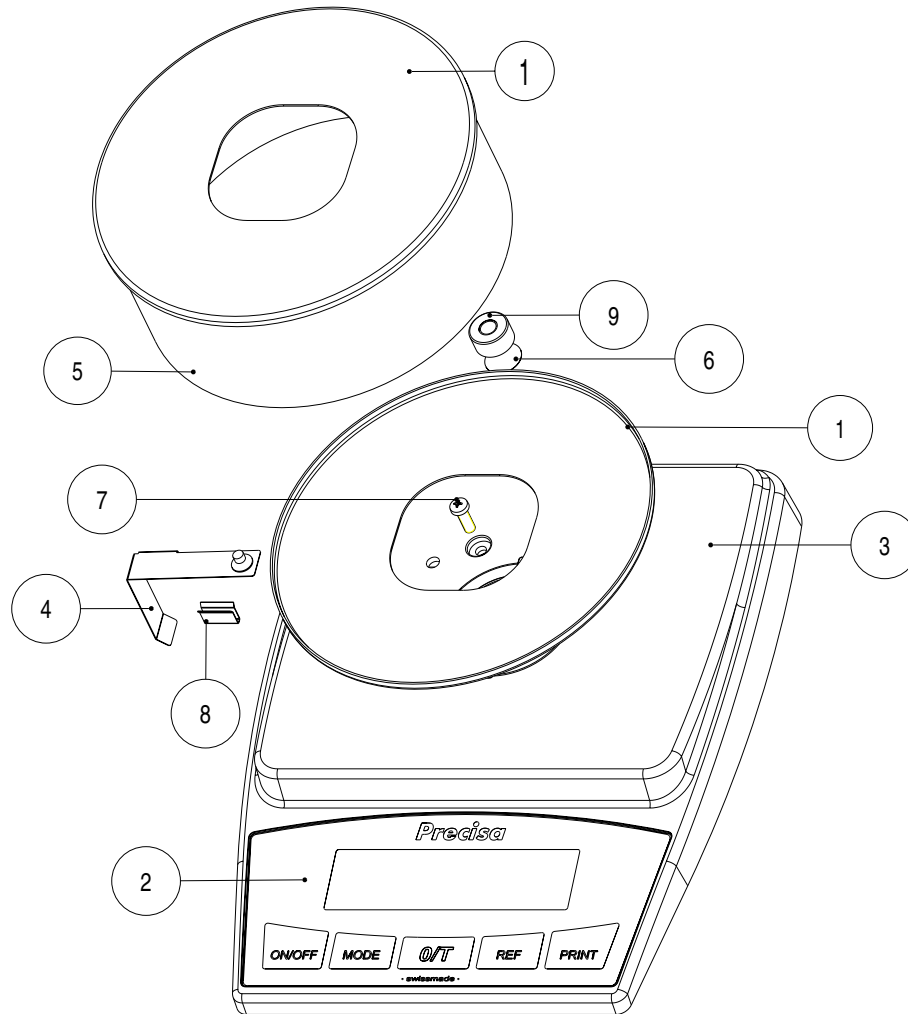
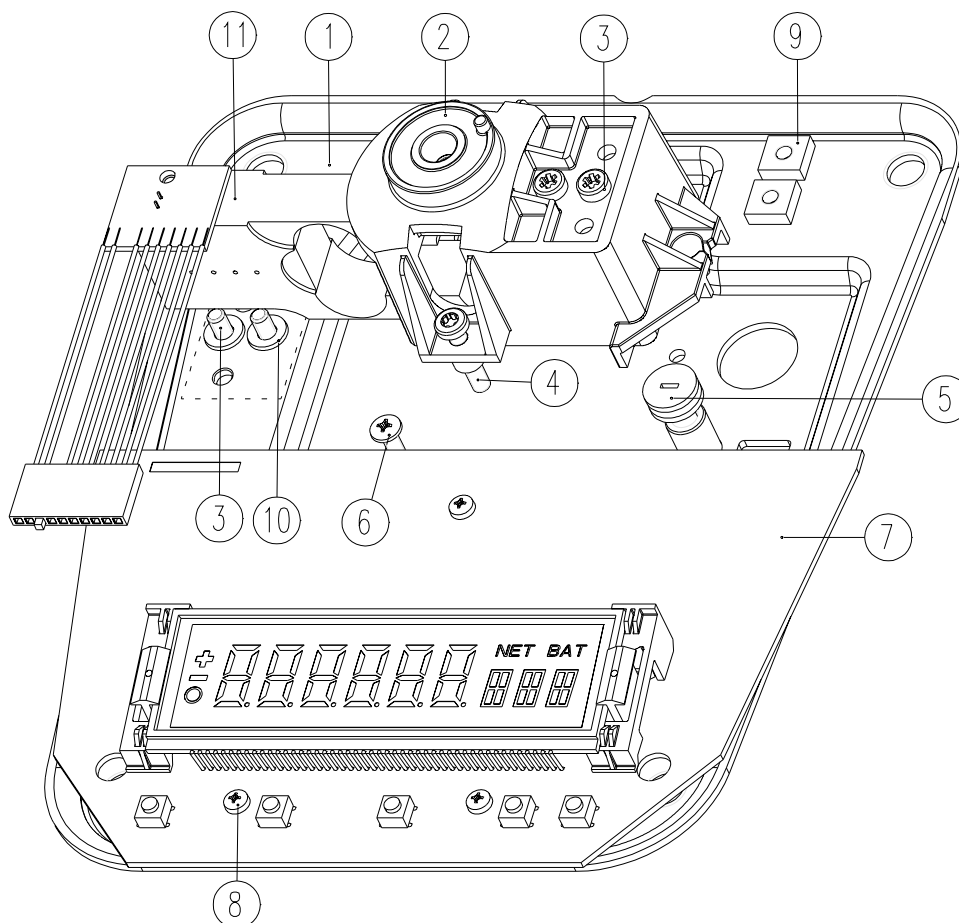
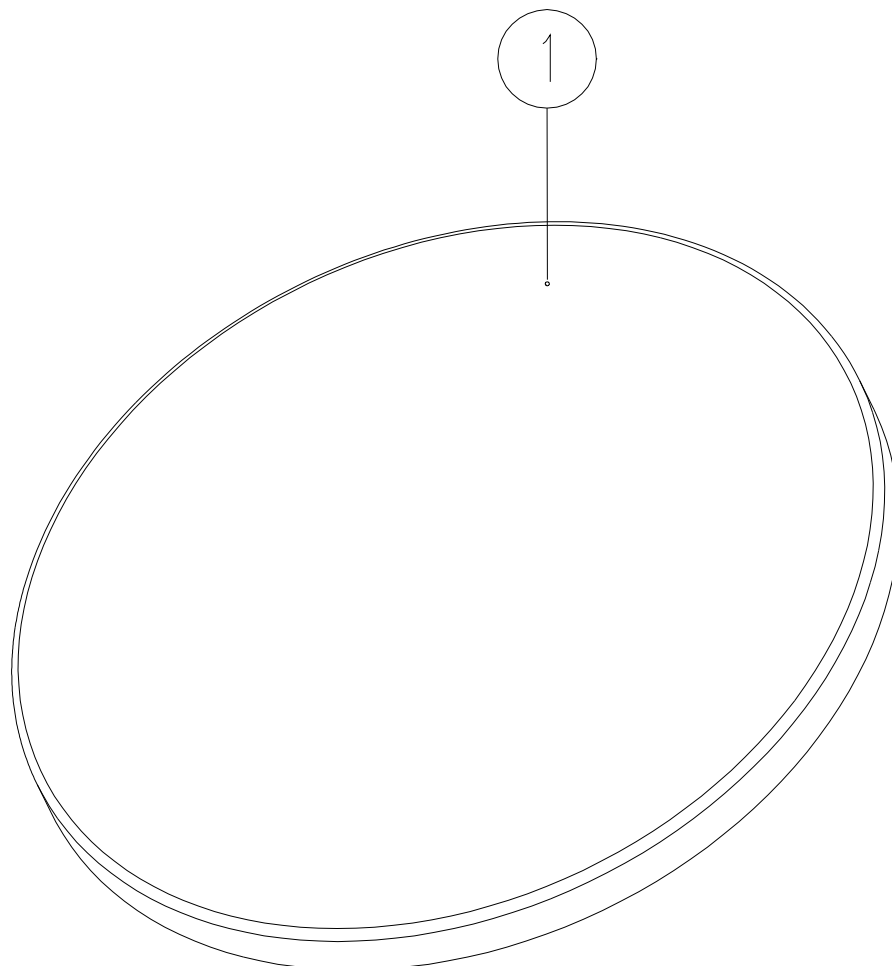
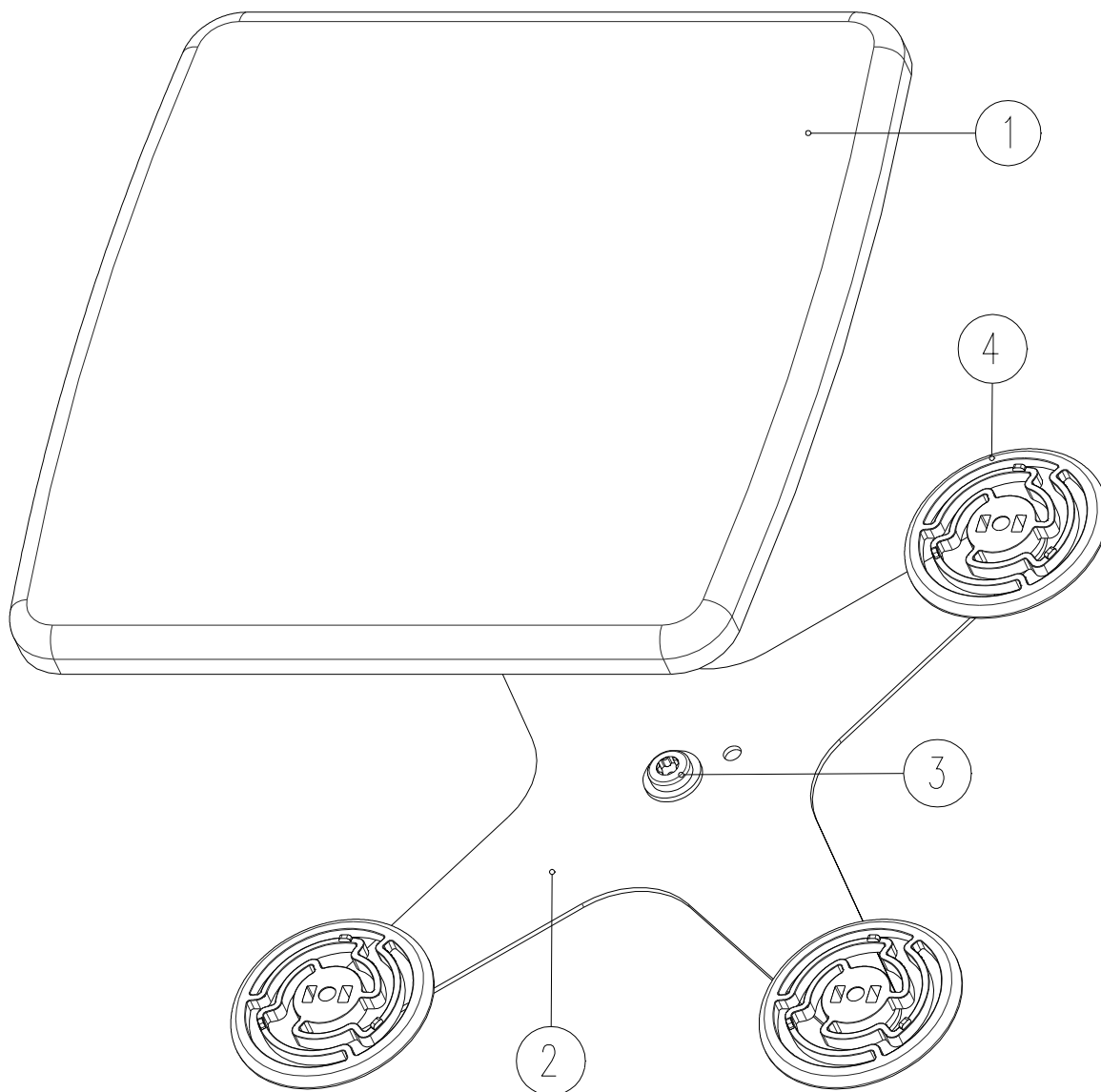
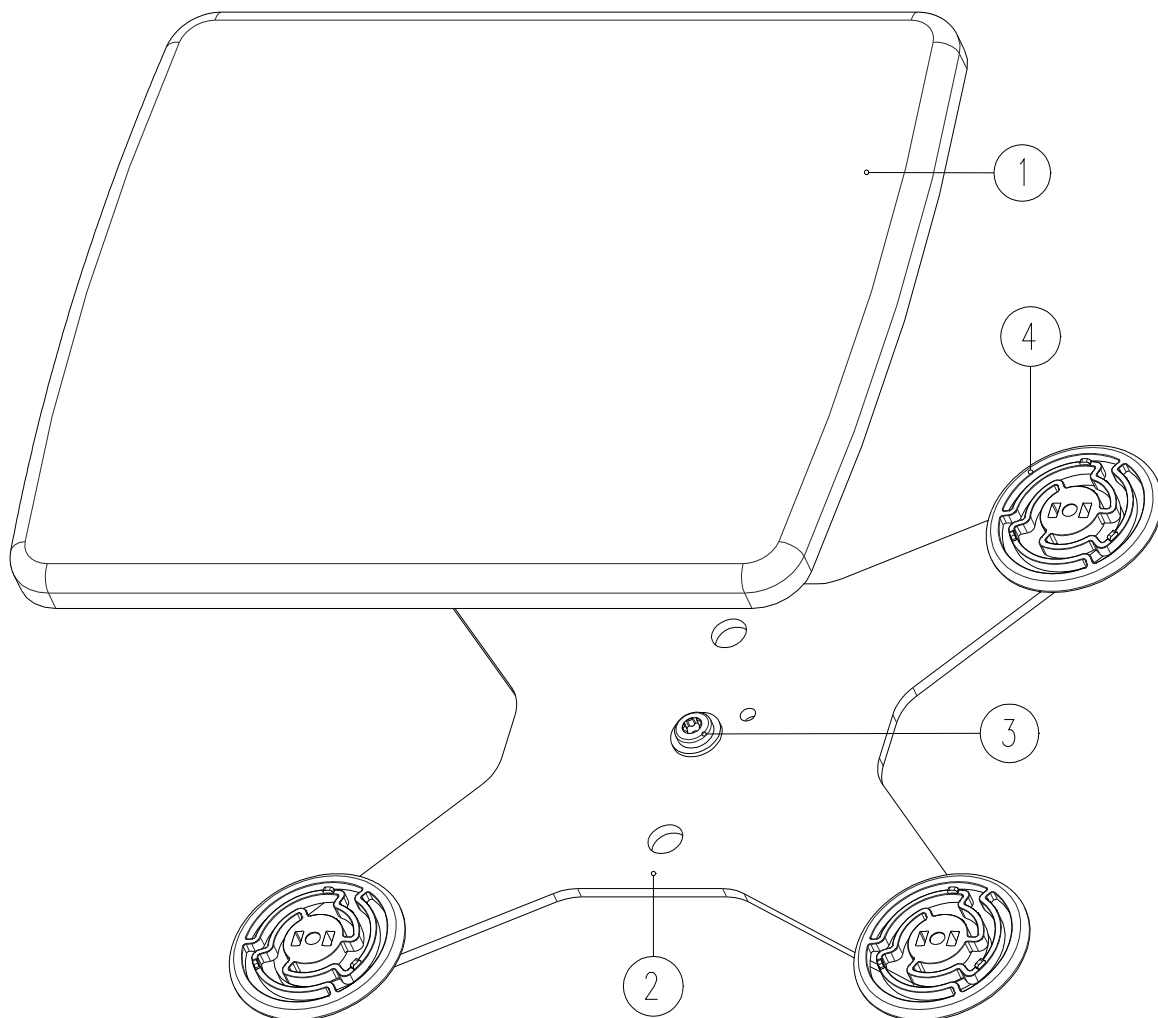
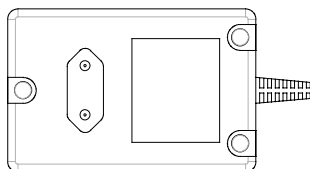
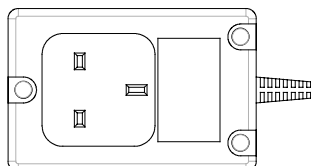
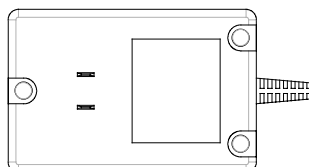
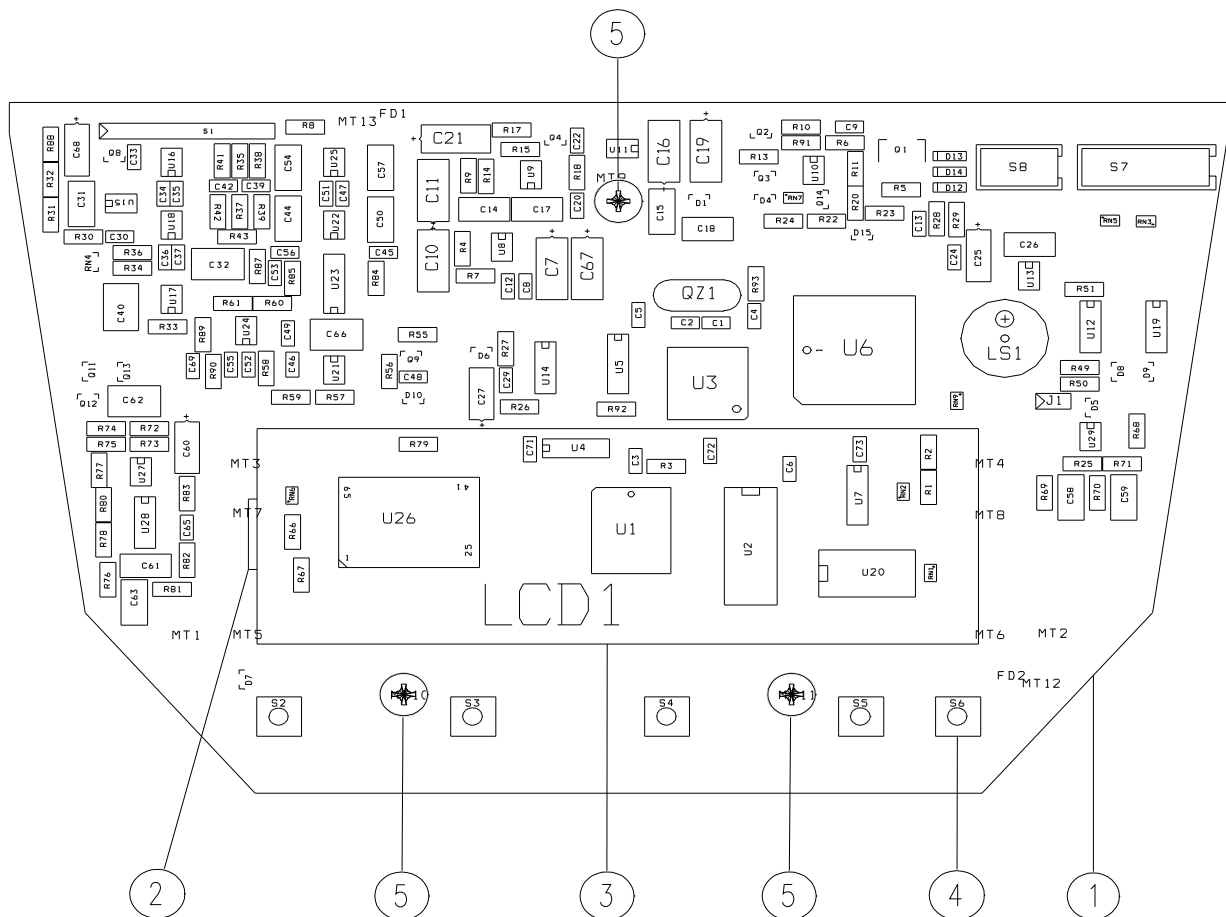
2. Main board

Important!

When the main board is replaced, there is the need to programme the flash with the **balance specific file „pee“**. The balance specific file „pee“ is an individual and specific data file and consequently only applicable for the concerned balance. For this reason you need to demand the file at Precisa Gravimetrics AG and download it into the flash-memory.

1. Open the balance, see page E2, chapter 27..
2. Remove the ribbon on the main print side, which leads from the connector print to the main print.
3. Remove the three screws (5) from the main print.
4. Lift up the main print.
5. Install the main print in reverse order.
6. Programme the flash-memory with software nnnn-nnn.PEE, see page E11.

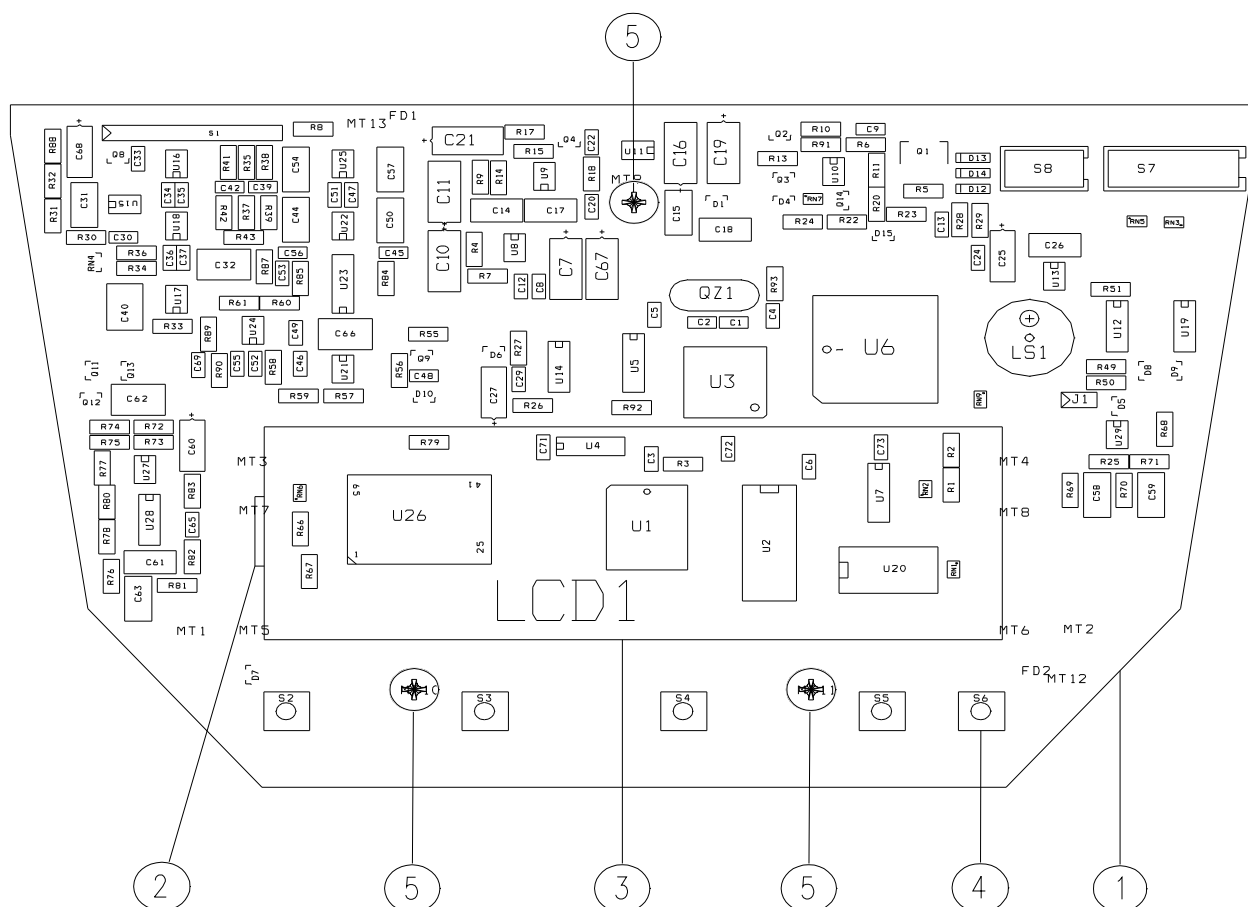Key-No.	Part-No.	Description
1	160-7205-020	Main board without LCD display
1	160-7200-020	Main board with LCD display
2	160-4008	Holder, for LCD display
3	PN 3260-021	LCD display
4	PN 3602-011	Micro switch

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Important!

By ordering, pay attention if the connector print version is equipped **without** interface for second display (Illustration 1) or **with** interface for second display (illustration 2).

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
Section E: Adjustment


1. Open up the balance

1. Disconnect the balance from the mains.
2. Take off the weighing pan.
3. Remove screw (1) from the pan holder/support by using the Line-head screwdriver included in the service-tool from the 320 series balance.
4. Use a No. 2 Philips screwdriver, to detach screw (2) at the top of the housing.
5. Use the same Philips screwdriver, to detach screws (3) at the bottom of the housing.
6. Lift up the housing's top.

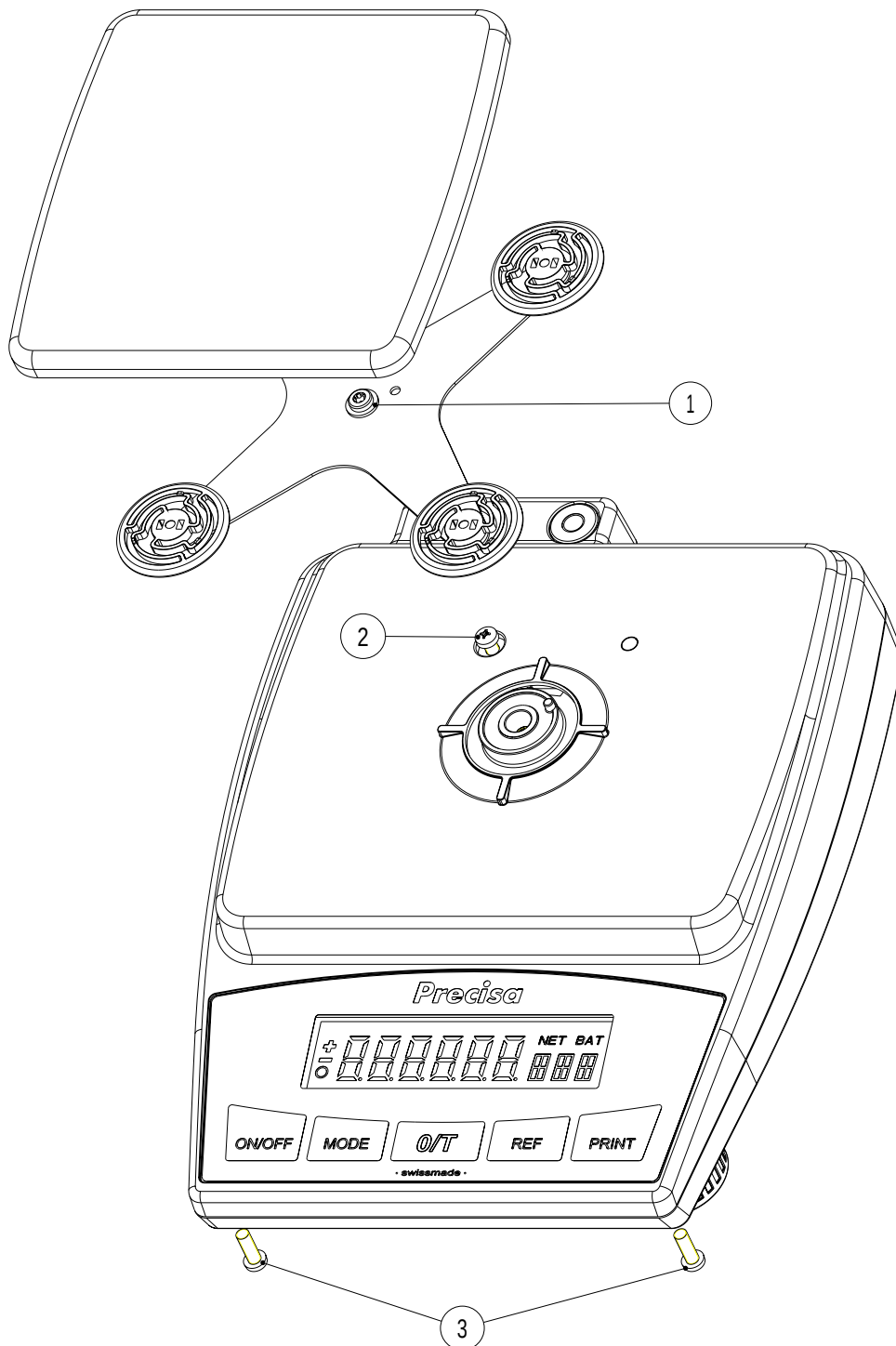2. Removal of the base plate

1. Open up the balance, see page E2.
2. Remove the ribbon cable (1) of connector board (3) and main board.
3. Remove screws (2) from connector board (3). **Pay attention to the battery clip!**
4. „M“ and „C“ balances are only equipped with levelling screws in front. „D“ and „G“ balances are also equipped with levelling screws at the rear.
-> Remove all levelling feet (5), thread sleeves (4) as well the fastening sleeves (6) with a hex socket wrench or an other similar tool.
5. Remove screw (7) from the base plate.
6. The base plate can now be lifted out.

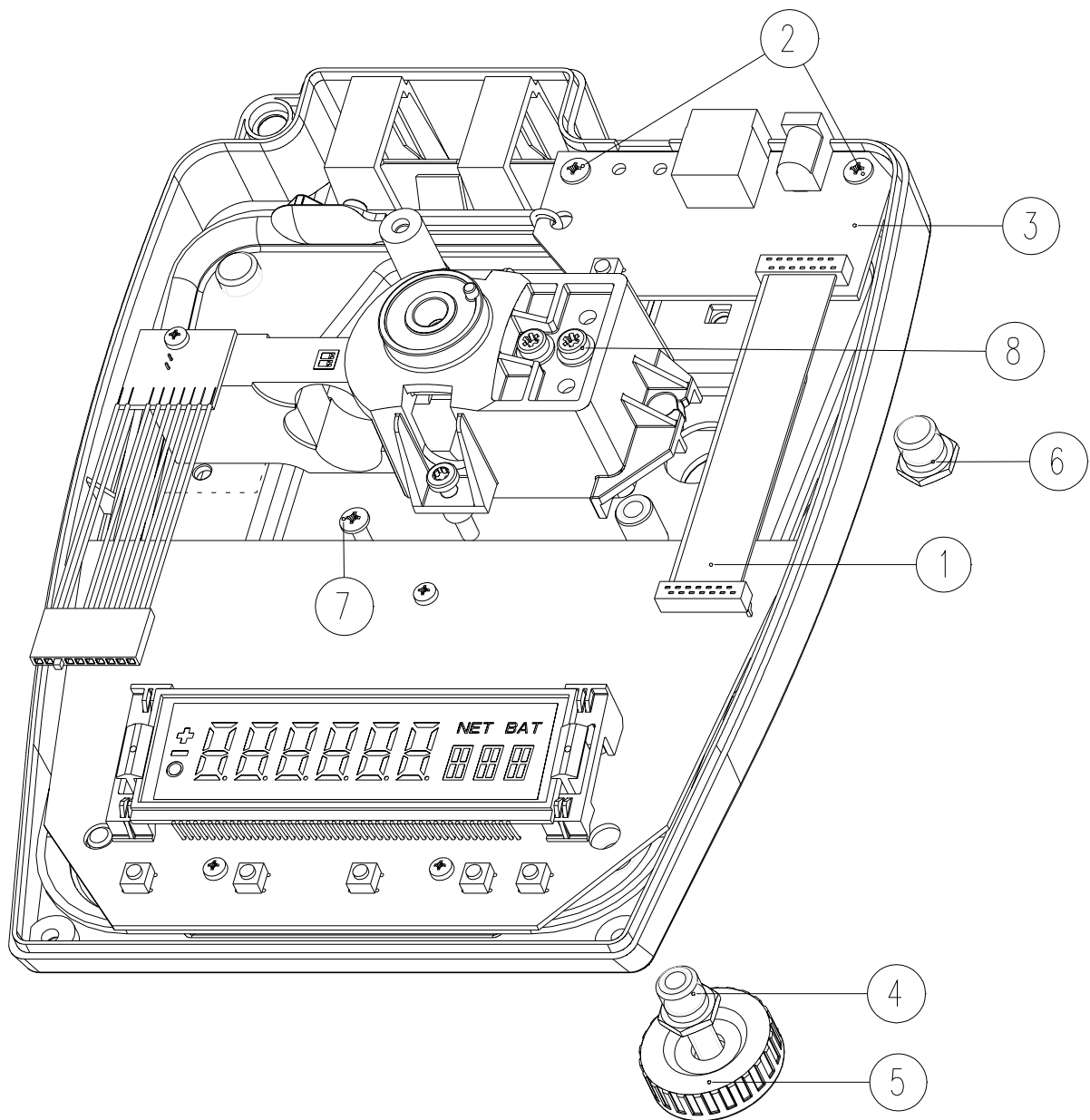3. Replacing the beam

3.1 Remove the beam


1. Open up the balance, see page E2.
2. Remove the base plate, see page E3.
3. Remove connecting cable (1). If the balance is equipped with accumulator or battery, also remove the ribbon cable connecting accumulator and battery with the main board.
4. Remove the underload adjusting screw (6).
5. **Important!** Do not remove the two overload adjusting screw, see page E5, leave them in the original manufacturer setting.
6. Remove the Line-Head screws which fasten the beam (3) on the underside of the base plate, with the Line-Head screw driver from the 320 series tool kit.
7. Lift up the beam from the base plate.
8. Two or three Line-Head screws (4), depending on the model, must be removed from the support piece (5).

3.2 Installing the beam

1. Align the beam parallel to the margin of the rear base plate edge.
2. Re-assemble the new beam in reverse order.
3. After replacement of the beam it is necessary to download the balance specific file „pee“ into the flash-memory, see page E11.


3.3 Adjustment and check of over- and underload

4. Put the balance in service mode, see page E5, chapter 30..
5. Adjustment and check of over- and underload, see page E6, chapter 31...

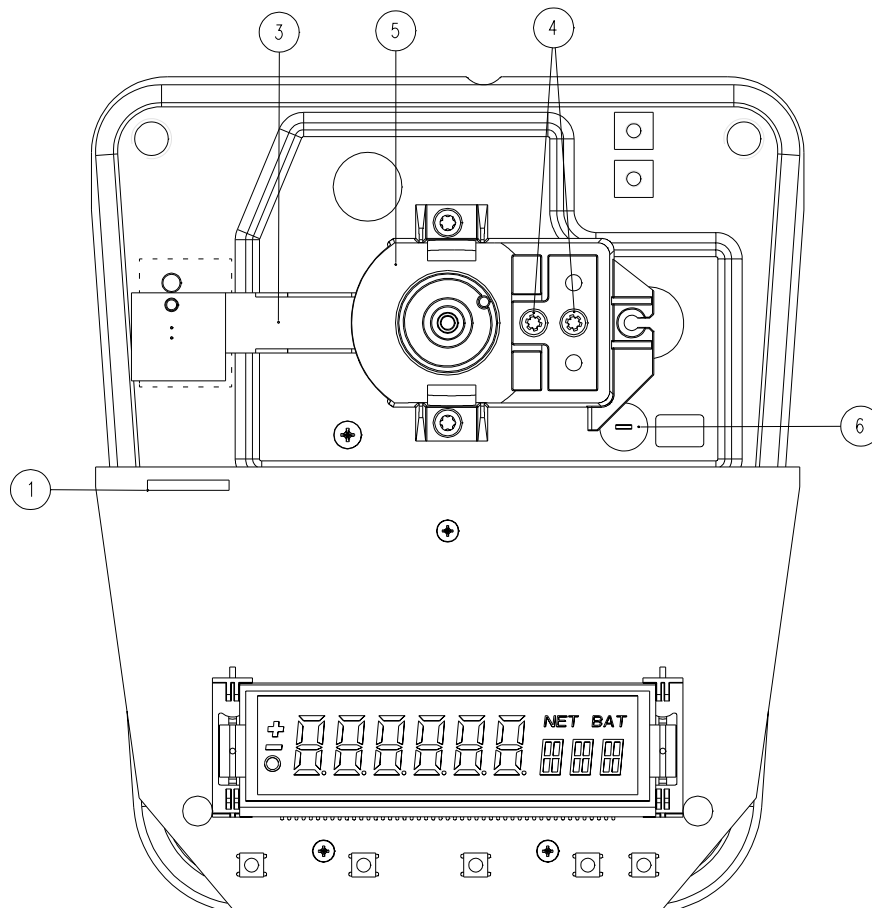4. Inspection / Adjustment

1. The top part of housing must be removed.
2. Disconnect the balance from the mains.
3. Remove the jumper.
4. Assemble the pan holder and fit the weighing pan.
5. Connect the balance to the mains. If the balance is equipped with accumulator or battery the necessary adjustments may be made without connecting the balance to the mains.
6. Press „ON/OFF“ , after the balance is switched on, the balance will be in service mode.

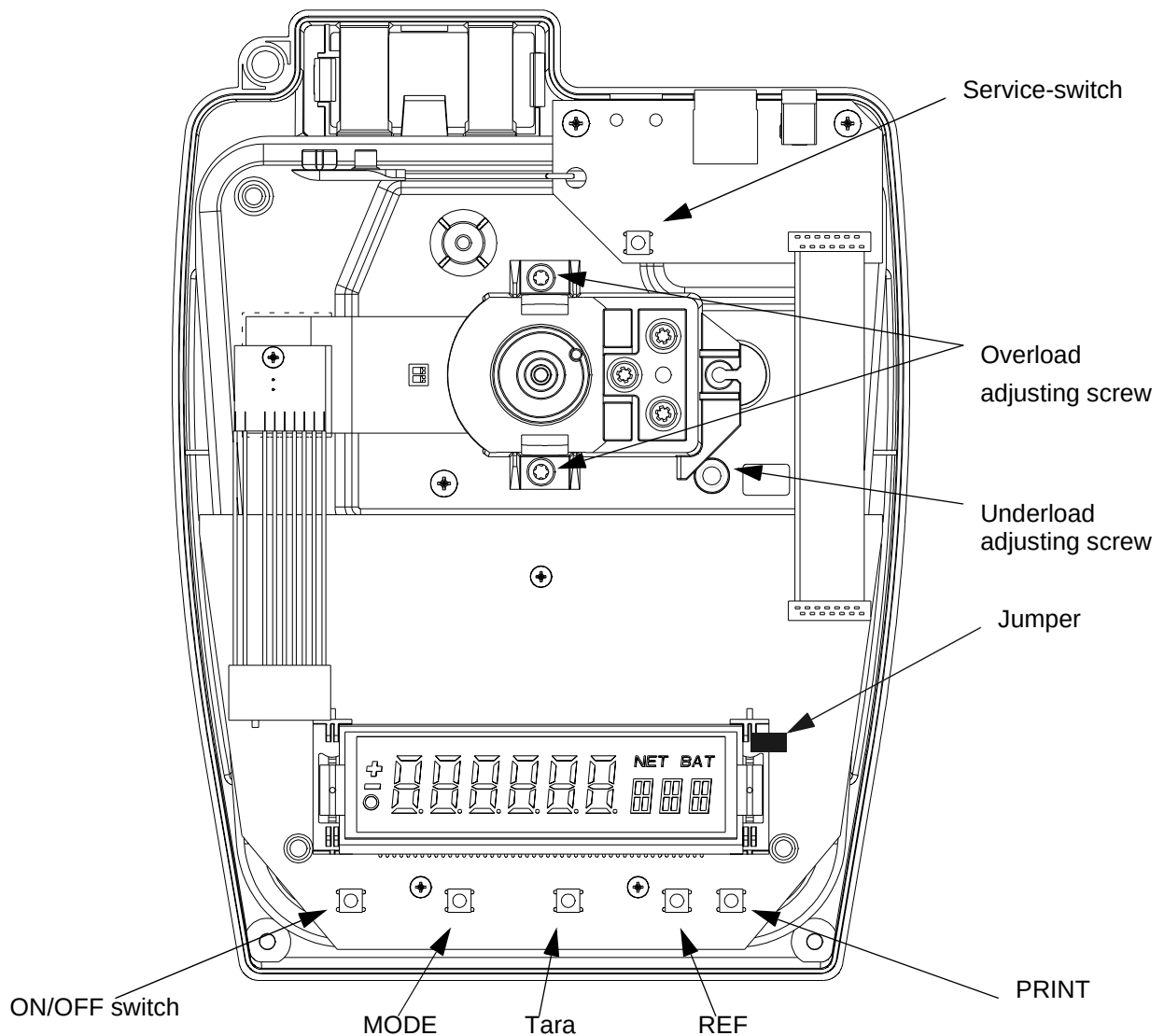5. Check and adjustment of over-, and underload protection

Balance must be in service mode, see page E5, chapter 30..

5.1 Checking overload

1. Determine overload check weight, see data sheet page E10.
2. Hold down the key „Mode“ until „TEST 1“ appears on the display.
3. Load the balance with the check weight. You may also press with your hand on the weighing pan.
-> With the check weight on the pan or the pressure you do exert with your hand on the pan, the tolerances for overload „from - till“ must be between the figures indicated in the data sheet page E10.
-> **Important!** If you exert the pressure with your hand **do not over press** the tolerance figures indicated in the data sheet. If you do over press, the beam may be damaged. Adjust the overload screw by turning clockwise till the screws touch the base plate and you can not any more overload over the tolerances indicated in the data sheet.
4. The adjustment of overload must be made with the two overload adjustment screws, see page E5, this is necessary to assure the screws are set to the same distance to the base plate. Turning clockwise overload will be reduced, turning anticlockwise overload increases.
5. **Important!** Take care that the adjustment is made with the **two** screws for overload. Turn only 1/4 in the direction you need to turn.
-> If the adjustment is not done with the two screws and both screws do not have the same distance to the base plate, then the beam may tip over and be damaged.

5.2 Checking under-load

1. The underload may be adjusted with the underload adjusting screw, see page E5.
2. The distance from the support piece to the stop screw needs to be adjusted to 0,5 mm with a feeler gauge.

6. Check corner load

The corner load adjustment on the spare beam has already been done on factory. So a later adjustment of corner load is on principle not necessary. Should an adjustment still be necessary, follow strictly the instruction on chapitre 32.1.

1. Determine the weight for corner load check, see page E10.
2. Put the weight in the middle of the pan. Point „M“, see illustration.
3. Take note of the displayed value.
4. Put the weight on points A, B, C, D, see illustration 1. Do not slide the weight to the alternative positions but instead lift it up to move it from point to point. If you have a round pan (M balance) you need to check the corner load as you do it with a square pan.

6.1 Adjustment of corner load


To be conscient of the effect of corner load change when filing off the beam do it only with light pressure on that way you will achieve the desired change in corner load.

1. Prepare the balance without the top of the housing but with the pan holder mounted on the support and the weighing pan on the pan holder.
2. Cover the main board with a sheet of paper or cloth to protect it from metal filings.

6.2 Example of corner load adjustment

1. Determine the check weight for corner load adjustment. See data sheet on page E10.
2. Put the weight on the centre of the pan, point „M“, see illustration 1.
3. Determine the lowest corner value between A, B, C, and D in relation to point „M“. The corner load adjustment starts at the highest negative deviation. on the corner with the major negative deviation.
4. Apply the file carefully on the relevant point A, B, C or D and file in the push direction only. See illustration 2 page E8.
5. Check the corner again as explained on chapter 32.

Illustration 1


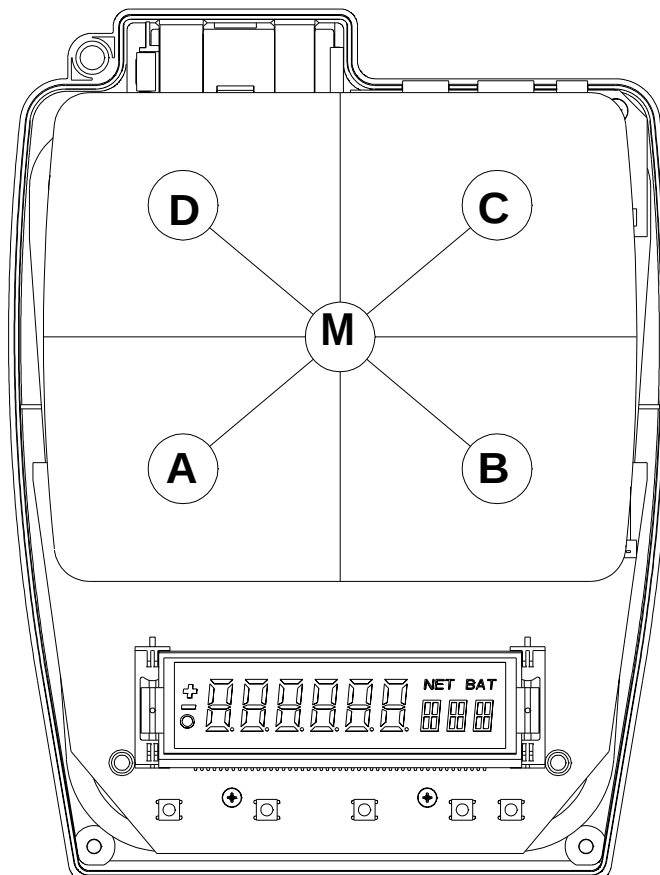
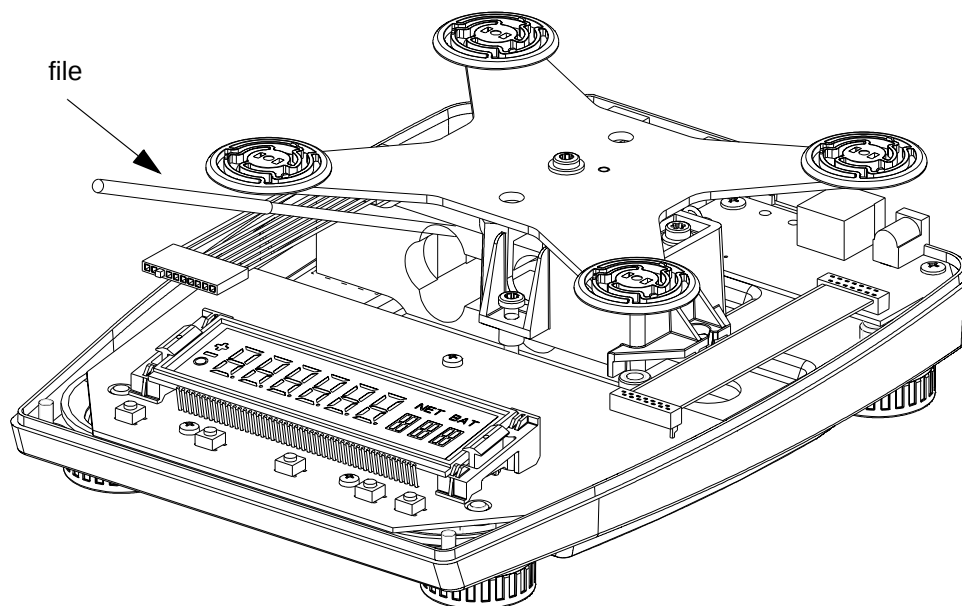
Illustration 2

front view

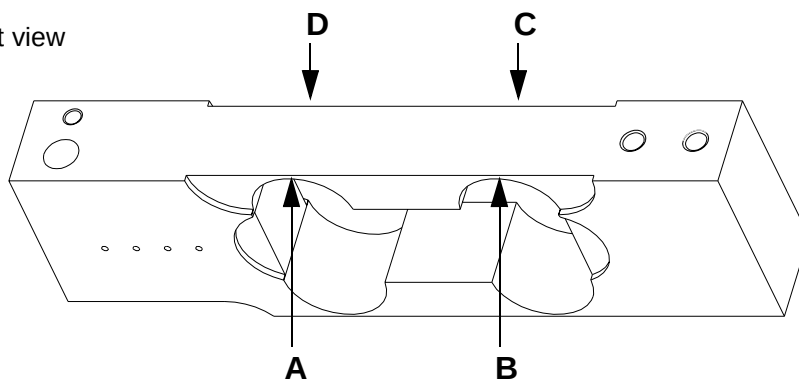rear view

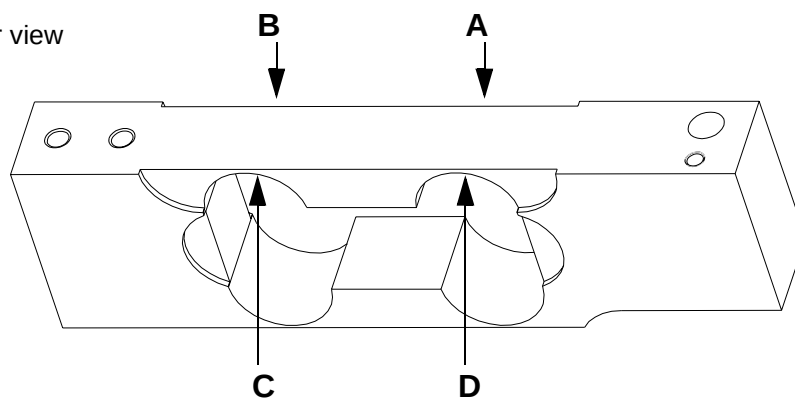6.3 Linearity of the balance

The balance must be in the service mode and have reached its correct operating temperature (this is achieved 15 min. after turning the balance on).

1. Use half and full load weights, see data sheet page E10.
2. Press and hold the „MODE“ key until „linear“ for linearity is displayed.
3. Assure that weighing pan is empty. Press the „REF“ key and hold it, till „0“ is displayed.
 - > release „REF“ key.
 - > the balance starts to establish the zero point and flashes. The zero point determination is finished as soon the flashing stops.
 - > The figures displayed are internal balance values and not relevant to this operation.
4. Put the half full load weight on the centre of the pan.
5. Press and hold the „REF“ key till the half weight specification appears.
 - > release „REF“ key.
 - > the balance starts to establish the half load weight and flashes. The half load is determined when the display stops flashing.
 - > The figures displayed are internal balance values and not relevant to this operation.
6. Put the full load weight on the centre of the pan.
7. Press and hold the „REF“ key till the full weight specification appears.
 - > release „REF“ key.
 - > the balance starts to establish the full load weight and flashes. The full load is determined when the display stops flashing.
 - > The figures displayed are internal balance values and not relevant to this operation.
 - > Full load may remain on the weighing pan.
8. Press and hold the „REF“ key till „SET“ appears on display.
 - > release the „REF“ key.
 - > The balance calculates the linearity - display is flashes - calculation is achieved and the operation is completed as soon as the full load capacity is indicated on the display.
9. Disconnect the balance from the mains.
10. Replace the jumper.
11. Re-assemble the balance again in reverse order.
12. Align the pan holder on the support when you fix the pan holder with the screw.

6.4 Final check of the adjustment

1. Switch on the balance.
2. The balance must have reached its correct temperature (this is achieved 15 min. after turning the balance on).
3. Level the balance with the levelling screws.
4. Calibrate the balance with the „T“ key.
5. Check corner load. See data sheet.
6. Check linearity. See data sheet.
7. Check reproducibility.

7. Programming the Flash-memory with the balance specific Software nnnn-nnn.PEE

Pre information:

To program the flash-memory with the balance specific software you need to ask at Precisa Gravimetrics AG for the **balance specific file „PEE“**, do that through the e-mail address service@precisa you will then get the file by e-mail. Each balance possess an individual and specific file, therefore you need to indicate the series number of the concerned balance. You find the series number on the rear of the balance on the white small sticker near the bubble level or inside the balance on the base plate rear on the left side.

7.1 Installing the EEPROM writer

1. Load the Precisa Gravimetrics AG homepage www.precisa.com
2. Select the field „Downloads“.
3. Select the field „Agents“ (Restricted Area).
4. Feed user name and codeword into the computer. Both access codes has been given by sales dpt. with flash no. 500.
5. Select the file „EEProm writer/“.
6. Click on the [EEWriterVnn.exe](#).
 - >download [EEWriterVnn.exe](#) on your PC/Laptop.
 - >define the folder where you want to install the [EEWriterVnn.exe](#).
 - >select the folder where you saved the [EEWriterVnn.exe](#).
 - >double click on the file to start the „unzip“ procedure.
 - >three files will be opened.
 - >click on the file „setup.exe“. You will be asked to install the download program [EEWriterVnn.exe](#).
 - >click on the big square with a PC icon.
 - >the „program group“ Precisa will be suggested, confirm with „Continue“. The successfully „set up“ will be confirmed with „Precisa EEPROM Writer Setup was completed successfully“.
7. select the windows bottom „start“
 - >select „Program“
 - >select „Precisa“
 - >select „[EEWriterVnn.exe](#)“.
 - >the EEPROM Writer will be loaded.

7.2 Programming the Flash-memory

1. Load the balance specific file sent by Precisa Gravimetrics AG on your PC/Laptop.
2. Balance is disconnected from the mains.
3. The housing's top is not mounted.
4. Remove jumper.
5. Connect the balance to the mains and switch on the balance pressing the „ON/OFF“ key.
If the balance is equipped with a battery or an accu then there is no need to connect the balance to the mains.
6. Connect the balance over the RS232 interface with your PC/Laptop's interface. The data cable you need for this connection is the same you've got with the tool kit for the 320 series balance.
7. Select on the „EEPromWriter“ tool, interface. Check if the right „port“ referred to your PC/Laptop has been selected. Change if necessary and confirm with „OK“.
8. Select „file“ and then „load“. Load the balance specific file supplied by Precisa Gravimetrics into the EEPROM writer.
9. Start the download by clicking on the green arrow.
10. The new programming of the flash-memory is achieved as soon you get the message „download successfully“ on the EEPROM writer tool.
11. Fix the jumper and the connection for battery and accu if balance was equipped with and assemble the balance in reverse order.
12. Do now the linearisation as described on page E9.
13. The balance is now again ready for operation.

8. Data sheet - adjustment and final control

Balance type	Weighing range (g)	Readability (g)	Reproducibility (g +/-)	Check weight (g) overload adjustment	Underload (g)	
					from	till
60M	60.000	0.001	0.002	150	90	140
100M	100.000	0.001	0.002	150	120	140
210C	210.00	0.01	0.01	1000	400	900
410C	410.00	0.01	0.01	1000	600	900
610C	610.00	0.01	0.01	1000	800	900
1000C	1020.00	0.01	0.02	1500	1100	1200
1100D	1100.0	0.1	0.1	10000	3000	9000
2100D	2100.0	0.1	0.1	10000	4000	9000
4100D	4100.0	0.1	0.1	10000	6000	9000
6100D	6100.0	0.1	0.1	10000	8000	9000
8100D	8100.0	0.1	0.2	10000	9000	10000
4100G	4100	1	1	20000	6000	9000
6100G	6100	1	1	20000	8000	9000
12100G	12100	1	1	20000	15000	18000

Balance type	Weights for linear. (g)	Tolerance linear. (g +/-)	Corner load weight (g)	Corner load tolerance (g +/-)	Recommended calibration weight
60M	2x 30	0.002	20	0.002	50
100M	2x 50	0.002	40	0.002	100
210C	2x 100	0.02	100	0.02	200
410C	2x 200	0.02	100	0.02	200
610C	2x 300	0.02	200	0.02	500
1000C	2x 500	0.02	300	0.02	1000
1100D	2x 5000	0.2	200	0.2	1000
2100D	2x 1000	0.2	1000	0.2	2000
4100D	2x 2000	0.2	1000	0.2	2000
6100D	2x 3000	0.2	2000	0.2	5000
8100D	2x 4000	0.2	2500	0.2	5000
4100G	2x 2000	2.0	1000	2.0	2000
6100G	2x 3000	2.0	2000	2.0	5000
12100G	2x 6000	2.0	5000	2.0	5000

9. Programming the Flash-memory with the balance program software snn-nnnn.BIN

9.1 Installing the Precisa „DownloadVnnn.exe“ program

1. Load Precisa Gravimetrics AG homepage www.precisa.com.
2. Select the field „Public Download Area“.
3. Select „software“.
4. Select the field „Tool for downloading“.
5. Click on [DownloadVnnn.exe](#).
 - >download [DownloadVnnn.exe](#) on your PC/Laptop.
 - >define the folder where you want to install the file [DownloadVnn.exe](#).
 - >select the folder where you like to store the program [DownloadVnn.exe](#).
 - >double click on the file to start the „unzip“ procedure.
 - >three files will open
 - >click on the file „setup.exe“. You will be asked to install the download program [DownloadVnn.exe](#).
 - >click on the big square with a PC icon.
 - >the „program group“ „Precisa“ will be suggested, confirm with „Continue“. The successfully „set up“ will be confirmed with „Precisa download V2.03 setup was completed successfully“
6. select the windows bottom „start“
 - >select „Program“.
 - >select „Precisa“.
 - >select „Precisa Download V2.03“.
 - >the program V2.03 is installing,ile you selected before and double click.

9.2 Programing the Flash-Memory

1. Load the Precisa Gravimetrics AG homepage www.precisa.com
2. Select „Public download area“.
3. Select the file [Balance BJ](#).
4. Select the software [snn-nnnn.BIN](#) and save it on your PC/Laptop.
5. Balance is disconnected from the mains.
6. Re-assemble the balance in reverse order.
7. Connect the blance to the mains and switch on the balance by pressing on „ON/OFF“ key. If the balance is equipped with a battery or an accu then ther is no need to connect the balance to the mains.
8. Connect the balance over the RS232 interface with your PC/Laptop's interface.
9. Connect the balance over the RS232 interface with your PC/Laptop's interface. The data cable which ist needed for this connection is the same in the tool kit for the 320 series balance.
10. Select on the [Download Vn.nn](#) program, interface.
11. Select „Auto setting“ and confirm with „OK“.
12. Select „file“ and then „load“. Load the balance program software snn-nnnn.BIN into the [Download Vn.nn](#) tool.
13. Start the download by clicking on the green arrow.
14. New programming of the flash-memory is achieved as soon as „download successfully“ appears on the display.
15. Do now the linearisation as described on page E9.
16. The balance is now again ready for operation.

10. Error messages

<i>ERROR-No.</i>	<i>ERROR responsibility</i>	<i>ERROR repair</i>
-------------------------	------------------------------------	----------------------------

10.1 Operating error (Error disappears with correct use or application)

ERROR 1	invalid program	load right program
ERROR 3	weight too small	use higher weight
ERROR 4	wrong anti-theft-code	use right code
ERROR 5	wrong anti-theft-code	use right code
ERROR 8	zero point outside range	
ERROR 9 / 0	calibration factor bad	repeat calibration
ERROR 9 / 1	zero point difference too high	repeat calibration
ERROR 9 / 2	difference of internal weight too high	repeat calibration
ERROR 9 / 3	calibration factor bad	repeat calibration
ERROR 30	no internal reference weight	acknowledge with TARE-key
ERROR 31	no external reference weight	acknowledge with TARE-key
ERROR 47	too much values	use less values
ERROR 48	other unit then first value	use the original unit
ERROR 49	too less values	use more values
ERROR 53	values out of range +/- 50%	use right weight
ERROR 60	minimum transformer area fallen below	
ERROR 61	maximum transformer area exceeded	

10.2 Fatal error (Balance stays still, error repair in service mode)

ERROR 10	no linearisation values	repeat linearisation
ERROR 11	no calibration factor	repeat calibration
ERROR 21	no SCS temperature compensation values	
ERROR 22	no temperature compensation values	
ERROR 23	no SCS temperature no-linearity compensation values	

10.3 Hardware error (Balance stays still)

ERROR 14	wrong model code	call service
ERROR 16	internal RAM destroyed	replace CPU
ERROR 17	external RAM destroyed	replace RAM
ERROR 18 / 1	wrong program check sum	load right program, replace FLASH
ERROR 18 / 2	no approved program build in approved balance	load right program
ERROR 18 / 8	FLASH empty	load right program, replace FLASH
ERROR 18 / 9	wrong controller check sum	replace CPU
ERROR 19	wrong code figure	load right code figure
ERROR 20	temperature measurement defective	call service

11. Discussion with the manufacturer

In discussions with the manufacturer's works the following information should be provided:

- **Balance type** e.g. XT 220A
- **Order number** e.g. 320-9201 (Label on the rear of the balance)
- **Sales number** e.g. N 43210 (Label on the rear of the balance)
- **Series number** e.g. 2400-100 (Label on the rear of the balance)
- **Program number** e.g. 351,00,00 (appears on the display when switching on the balance)
- **As precise a description of the defect as possible**
- **Weighing results**

PRECISA Gravimetrics AG
Moosmattstrasse 32
Po-Box 352
CH-8953 Dietikon
Tel. +41 1 744 28 28
Fax. +41 1 744 28 38
e-mail: service@precisa.ch
Internet: <http://www.precisa.com>

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Hauptprint 2
Hersteller 14
Housing's bottom 2
Housing's top 3, 4

M

Main 3
manufacturer 14
Mechanics 1
Mechanik 1

N

Netzadapter 12

O

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P

Plate 7
Power 12

Power supply 12

R

Rücksprache 14

T

Top 3, 4

W

Waagschale 135x135 10
Waagschale 140x170 11
Weighing 9, 10, 11
Weighing pan 135x135 10
Weighing pan 140x170 9, 11

Z

Zubehör 1