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SHOULDER MASSAGER

IMPORTANT SAFETY PRECAUTIONS

- The following persons must consult their physician before using the product.
 - If you have a pacemaker or similar implanted medical device.
 - If you have a heart condition.
 - If you have a strained back.
 - If you are undergoing medical or paramedical treatment.
 - If you have weak or injured bones, especially back bones.
 - During pregnancy or shortly after delivery
 - If you have osteoporosis.
 - If you have blood circulation problems or thrombosis.
 - If you have a skin infection.
 - After surgery.
 - If you have a physical abnormality.
 - If you have a spinal disorder or injury caused by a past accident or illness.
- Do not let water spill on or get into the product.
- Water and electricity is a dangerous combination. Do not use the product in wet surroundings (e.g. in the bathroom or near a shower or swimming pool).
- Use the product only at the designed voltage.
- If the power plug, main cord or the product is damaged, please stop using the product.
- To prevent electricity hazard, do not open the product's outer cover.
- This product is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Always return the appliance to our service centre for examination or repair if the product or its power cord is broken for safety reasons. Do not attempt to repair the appliance yourself, otherwise your guarantee becomes invalid.
- This appliance is not intended for use by children. Children should be supervised to ensure that they do not play with the appliance.
- Do not use the appliance if the outer cover is damaged or in case of mechanical or electronic failure. Contact our service centre immediately.
- Do not use the appliance in a room with high humidity and high temperature or next to a hot object (e.g. Heater).
- Do not put medical or cosmetic creams or lotions on the appliance or on your back before or during use of the appliance.
- Do not shake or drop the appliance.
- Do not damage, bend, pull, twist or bundle the power cord. Protect it from being placed under heavy objects or caught between objects while moving the product. A damaged cord may lead to fire or electric shock.
- The product should be started up and operated in the presence of people. Please pull up the power cord when not using it for a long time or with the absence of people.
- The continuous operating time of the machine should be no more than 30 minutes; restart the machine after a rest of 20 minutes.
- In case discomfort or abnormality occurs during application, please stop using it immediately.
- The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.

SUITABLE FOR

- Computer users; · Office workers; · Housewives; · Athletes;
- People with shoulder or back pain due to long-time use of arms;
- People with muscle fatigue and pull due to excessive exercise;
- Elderly people with unsmooth blood circulation in the back.

BENEFITS

- Relieve shoulder pain due to stiff neck during sleep.
- Relieve pain due to a long unsmooth blood circulation.
- Relieve fatigue due to excessive exercise.


- Relax the tense muscle.
- Helps to relieve stress.

FOR OPTIMUM RESULTS

- Use the product 2~3 times per day with 5~15 minutes per time to eliminate back pain;
- People operating computers or using arms for long time use it for 5~15 minutes to relax back muscles, ameliorate the blood circulation in the back and prevent injury of back muscle;
- The product is with temperature protection device. When the temperature exceeds safety limits, the machine will stop operating automatically. Please unplug the power supply for cooling before using it again.
- Please have a cup of warm water after using the product to enhance metabolism.

FUNCTION KEYS AND OPERATING INSTRUCTIONS

1. Function Keys

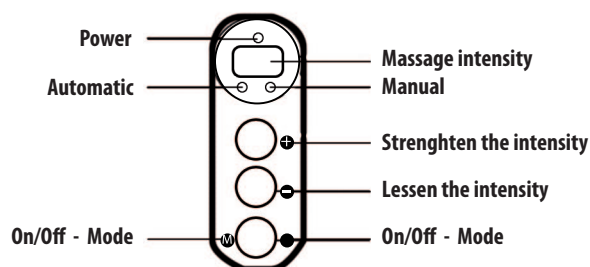- Press "ON/OFF/MODE" for two seconds to turn on
- When plugged, Power indicator flickers.

2. Operating Instructions

- Put the mains plug in the wall socket.
- Press "ON/OFF/MODE" for two seconds to turn on, the machine will enter the preset 6 massage modes, one auto mode and five manual modes. It starts operating from the first massage mode (auto mode) until the 6th mode in cycle. Press "ON/OFF/MODE" to switch the mode.
- This product is set with 6 level intensity, when you turn on, its preset intensity is the 3th level. You can adjust the massage intensity simply by press the "+" , "-" button. The LED indicates the current massage intensity.
- The massage session lasts 8 minutes. After the session, the appliance automatically switches to standby mode.
- If you want to stop the massage session, press the ON/OFF/MODE BUTTON for two seconds, the appliance will stop and enter into standby mode.
- After the massage session, remove the mains plug from the wall socket.

POSSIBLE USE OF MASSAGER

Shoulders - Upper back - Lower back - Hips - Thighs.

MAINTENANCE

Cleaning

- Never use scouring pads, abrasive cleaning agents or aggressive liquids such as petrol or acetone to clean the appliance.
- Do not immerse the appliance in water.
- Remove the mains plug from the wall socket.
- Clean the control panel and cords with a dry cloth.
- Do not use a moist cloth or a cleaning agent.

Storage

- Remove the mains plug from the wall socket.
- Store the appliance in a cool and dry place away from the sunlight.
- Do not put heavy objects on top of the appliance.

TROUBLESHOOTING

Condition	Possible Causes	Remedy
The power cannot be connected	Loosen power plug bad contact	Please replug the power plug
	You have used the appliance continuously for 25 minutes, causing the overheat protection to switch off the appliance	Unplug the appliance and let it cool down for 30 minutes before you see it again.
The product does not work	The product has just finished its 8 minutes massage and has automatically switched to standby mode	If you want to massage again, press the ON/OFF button.
	An excessive weight has been put on the product (e.g. someone has stood on the appliance or a heavy object has been put on it)	Unplug the power plug, remove the excessive weight if it is still on the product again and switch it on.
The product is warm.	Over temperature of the product	Unplug the power plug and let it cool down. Always unplug the power plug after use.
A sound comes out of the product	This sound is caused by the mechanical structure of the product	This sound is normal and does not mean that the product is malfunctioning

SPECIFICATION

Unit party : 760x480x65mm
 Net weight : 1.85Kg
 Power Consumption : 50W
 Power Input : AC230V ~ 50Hz
 International Compliance: CE FCC RoHS

WORKING ENVIRONMENT

Temperature : 0° ~ 40°
 Humidity: ≤85%
 Atmospheric pressure: 86kPa ~ 106kPa

WARRANTY

The warranty is valid for normal use of the device as defined in the user's guide. Excluded from this warranty, damage due to a cause unrelated to the device. In particular, the warranty does not apply if the unit has been damaged following a collision or a fall, a false move, a connection that does not comply with the instructions or requirements specified in this manual, the effect of lightning overvoltage sector, inadequate protection against heat, moisture or frost. This warranty gives you specific rights, and you may also have other rights that vary from state to state. If for any reason this product needs to be returned to the store where purchased it must be packed in its original carton/package. The manufacturer and distributor has the right to modify the apparatus (technically and aesthetically) and/or replace some parts, whereas necessary, to always offer a safe and reliable product.

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