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INN
904

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RECHARGEABLE BATTERY TOOTHBRUSH

USER MANUAL

The electric toothbrush PASTA DI CAPITANO is a device with sonic technology, it has an operating system of 30,000 vibrations per minute. The use of this toothbrush ensures a correct operation of brushing, to avoid disturbances due to excessive pressure or improper use of the bristles, fights the formation of bacterial plaque and improves dental hygiene and oral cavity.

Before starting using the appliance, please read carefully all the instructions below. It is advisable to save this instruction manual for a future use.

IMPORTANT SAFETY INSTRUCTION


When using electrical products, especially when children are present, the basic safety precautions below should always be followed.

READ ALL INSTRUCTIONS BEFORE USING:

- Before connecting the appliance, check the voltage shown on the rating label corresponds to the local mains voltage.
 - Do not use this product if the power cord or plug are damaged, if it is not working properly, if it has been dropped or damaged, or dropped into any liquid. In case of doubts, please address to a specialised centre.
 - If the charging station falls into water, do not try to reach it into the water, but pull out the power plug immediately instead!
 - Do not immerse the charger into water or any other liquid.
 - Always pull the power cord out of the outlet before cleaning the brush stand or the charging station.
 - Keep the unit away from heat, direct sunlight, moisture, sharp objects and similar ones.
 - Do not immerse the device in water or other liquid.
 - Do not handle with wet hands.
 - Do not use while bathing.
 - Do not place or store product where it can fall or be pulled into a tub, shower stall or sink.
 - Use only with the charging station and accessories supplied.
 - Do not insert any foreign object into any opening on the product.
 - Make sure that the power cord does not touch hot surfaces. Do not wrap the power cord around the device.
 - Always pulling from the plug and not from the cable.
 - Do not dismantle this product at any time other than when disposing of it.
 - Do not use outdoors or operate where aerosol products are being used or oxygen is being used or administered. This is a fire hazard.
 - Use of this product may interfere with the operation of pacemakers or other internal medical devices. If you use one of these devices you should consult your physician.
 - Use this product only as indicated in these instructions or as recommended by your dental professional.
 - Follow the instructions for use, all uses other than described in this manual are considered improper use. The manufacturer declines all responsibilities for damages to persons or properties resulting from improper use of the device.
- This appliance is intended for cleaning teeth and gums. Do not use it for purposes other than those intended. If you feel pain during use, contact your doctor.
- Keep the device out of reach of children.
 - The toothbrush is not suitable for use by children under 3 years of age.
 - This device is not a toy. This product is not intended to be used by children without direct adult supervision.
- CAUTION:** Do not leave the battery compartment open in the presence of children, choking hazard with the ingestion of the batteries.
- The appliance is not intended for use by persons (including children) with mental, physical or sensory capabilities reduced, except in cases in which such persons have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
 - Keep all packing material (plastic bags, cardboard, polystyrene, etc..) out of reach of the children, since they are potential sources of danger and must be disposed of as provided by law.
 - The product is intended for household use only, do not use the product for commercial purposes.

ATTENTION:

- Do not clean the head attachment, the handle, the charger station, the charger station cover or sanitiser in the dishwasher.
- If during the last 2 months you have had interventions on the teeth or gums, ask the advice of your dentist before using the brush.
- Consult your dentist if excessive bleeding after use of the device or if bleeding continues for more than a week.
- In case of oral problems, consult your doctor before using an electric toothbrush Pasta Del Capitano.
- Do not use the appliance with crushed or bent bristles. Replace the brush head every two months and when ever appear signs of wear.
- Use only the heads recommended by the manufacturer.
- If your toothpaste contains peroxide, baking soda or bicarbonate (common in whitening toothpastes), then roughly clean the print head with soap and water after each use, to prevent damage to the plastic.

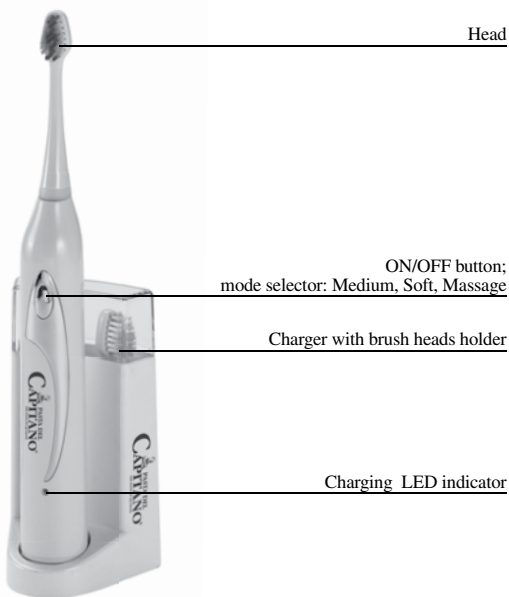PRODUCT DESCRIPTION**HOW TO CHARGE THE TOOTHBRUSH**

Charge the toothbrush for approximately 16 hours before the first use.

Insert the plug into the socket.

Place the switched off toothbrush on the charging station. The charging indicator lamp lights on and shows up the charging process.

A full recharge takes about 24 hours

Notes:

Remove the charger station power plug from the power outlet after the charging has completed.

Do not overcharge the toothbrush batteries)

A slight warming of the handle is normal during charging.

The charging takes place when the toothbrush is on the charging station and the toothbrush is switched off.

To extend the lifetime of the rechargeable batteries, it is recommend to drain them completely from time to time by using the toothbrush, and not to recharge them until the cleaning power is noticeably reduced.

To change the mode of use press several times the ON / OFF button.

To switch the device off press the ON/OFF button repeatedly until the device stops.

HOW TO INSERT THE HEAD ATTACHMENT

The brushes can be easily exchanged by inserting and removing them from the vibration shaft located in the upper part of the handle.

Pull up to remove the head from the handle. Insert a new head attachment until it clicks into place.

A gap must remain below the colored ring of the head for the free movement of the brush.

HOW TO USE THE BRUSH

- Moisten the bristles, apply toothpaste and switch the device on only after you have placed the brush in your mouth.
- Select the operating mode you want by pressing the ON/OFF button:
 - Daily cleaning - press the button once. This setting helps to clean in depth the spaces between the teeth and gums.
 - Gentle cleaning - press the button 2 times. This setting ensures a deep but gentle cleaning and is particularly suitable in the case of sensitive teeth.
 - Massage mode - press the button 3 times. This setting is particularly suitable for massaging and stimulating the gums.
 - Press the ON/OFF button to switch the device off.
- Swipe the brush slowly on each tooth, without pressing too hard toothbrush against the teeth. Clean the surface of the teeth, the inner sides and the edges of the teeth. It is recommended to brush teeth at least 1 minute.
- The device switches off automatically after 2 minutes (timer function)
- Do not remove the brush from your mouth until after the device has been switched off.

Note: In the first days of using an electric toothbrush some users may experience slight bleeding of the gums. This phenomenon should disappear in few days. If bleeding persists, please consult your doctor

To obtain a long-lasting, optimal cleaning effect, the brush attachments should be replaced every 2 months or whenever the bristles get damaged, stained or worn folded.

OBSERVANCE OF HYGIENE REQUIREMENTS

Thoroughly rinse the brush head after each use.

Remove the head from the handle and rinse under hot water, then dry before placing it on the handle. The head should be completely dry before you put it in the vibrating shaft.

Clean the handle with a damp cloth, then dry.

When necessary, wipe off the charging station with a moist cloth and rub it dry. Do not immerse in water! The station must be completely dry before you are permitted to connect it to the power line again.

Note: To keep the brush attachments in an hygienically perfect condition, it is recommended that they are put in a denture cleaning solution approximately evry 2-4 weeks in addition to the normal cleaning.

STORAGE

If the device is not used for a long time it is recommended to unplug it from the wall outlet and **store it in a dry place** away from heat sources.

Store the device well cleaned out of the reach of children. Do not wrap the power cord around the device.

TECHNICAL DATA:

Charging station:

Power: NiMH rechargeable batteries included (2x1.2V 3A 600mAh) - 230V~50Hz

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Handle:

Operating time per battery charged: approx. 70 minutes.

Charging time: first charging 16 hours; subsequent charging 12 hours.

Cleaning movements: 30.000/min., with a high amplitude (level 1); with a low amplitude (level 2); pulsing for massage (level 3).

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Read carefully the instruction manual.

This product complies with all the applicable EU directives.

This appliance is placed on the market after 13th August 2005 and must be separately collected.

Separate Collection- WEEE (Waste of Electric and Electronic equipment)

At the end of its use don't dispose the appliance, including removable parts and accessories, together with the other urban waste, but conformably to 2002/96/EC. Since it has to be treated apart from home waste, you must take it to a differential collection centre specific for electric and electronic equipment (EEE). Otherwise you can give it back to the retailer as you buy an equivalent apparatus. There will be severe sanctions in case of transgression. The correct disposal helps prevent potentially negative consequences on the environment and human health.

The batteries used in this device must be disposed of in the special bins at the end of their life. Please inform yourself about the local rules on separate collection of batteries. The correct disposal of batteries helps preventing potentially negative consequences on the environment and human health.

REMOVING THE RECHARGEABLE BATTERIES:

Before scrapping the device, you must remove the rechargeable batteries and dispose of them separately:

- Remove the brush.
- Loosen the screw on the base of the handle.
- Pry the base off.
- Turn the device so that it faces downwards, and place the vibration shaft on a hard surface.
- Pry off the two tabs to the casing middle that you can now see.
- Press the casing down.
- Now remove the rechargeable batteries from the compartment and use a pair of pliers to remove the cable that is soldered on.

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Prodotto sotto licenza ufficiale di **Farmaceutici Dott. Ciccarelli S.p.A.**