**IMPORTANTE**

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**IMPORTANT**

Before using the appliance read carefully all the instructions contained into this and all others documents attached, to know all the warnings, safety cautions, use instructions and illustrations. Remain this manual for future references.

Directive 2002/96/CE (waste electrical and electronic equipment – WEEE): information for users.

In accordance with European Directives 2002/95/CE, 2002/96/CE and 2003/108/CE on the restriction of use of dangerous substances in electric and electronic equipment as well as their waste disposal, this appliance is conformed.

The crossed-out wastebasket symbol on the appliance or on its packing means that at the end of its useful lifespan, the product must be collected separately from other waste. The user is responsible for delivering the appliance to appropriate collection facilities of electrical and electronic equipment waste, or for returning it to the retailer, buying a new equivalent appliance, on a one for one basis. Appropriate separate collection to permit recycling, treatment and environmentally compatible disposal helps to prevent negative impact on the environment and human health and promotes recycling of the material for making up the product. Abusive disposal of the product by the user involves application of the administrative sanctions according to the law in force.

**IMPORTANT**

Before using the appliance read carefully all the instructions contained into this and all others documents attached, to know all the warnings, safety cautions, use instructions and illustrations. Remain this manual for future references.

Directive 2002/96/CE (waste electrical and electronic equipment – WEEE): information for users.

In accordance with European Directives 2002/95/CE, 2002/96/CE and 2003/108/CE on the restriction of use of dangerous substances in electric and electronic equipment as well as their waste disposal, this appliance is conformed.

The crossed-out wastebasket symbol on the appliance or on its packing means that at the end of its useful lifespan, the product must be collected separately from other waste. The user is responsible for delivering the appliance to appropriate collection facilities of electrical and electronic equipment waste, or for returning it to the retailer, buying a new equivalent appliance, on a one for one basis. Appropriate separate collection to permit recycling, treatment and environmentally compatible disposal helps to prevent negative impact on the environment and human health and promotes recycling of the material for making up the product. Abusive disposal of the product by the user involves application of the administrative sanctions according to the law in force.

Information contained in this document have been carefully drew up and controlled. Nevertheless JOYCARE SPA Socio Unico does not accept any responsibility regarding possible mistakes. Information contained in this document, let alone the specifications of the referring products are subject to changes without notice. JOYCARE SPA Socio Unico does not accept any responsibility for the use or the application of the products described in

this document. If differently non specified, every refer to names, products, data, publishers and companies is purely accidental and it is only in order to illustrate the use of the JOYCARE SPA Socio Unico products. Names, places, marks quoted in case belong to the respective owners. All the rights are reserved and any parts of this document could be reproduced in any way and for any purpose without the wrote authorization of JOYCARE SPA Socio Unico .

### **Congratulations for the purpose of your PERFACTION!**

PERFACTION is an original and innovative device for the body care that synchronizes the vacuum method to your heartbeat, to maximize the action against the cellulite slight blemishes. Today thanks to the patented system, you can try out at your home an absolutely revolutionary treatment form that involves both the tissues superficial layer and the deeper, acting on the lymphatic system and stimulating its regular functions affected from the superfluous fat, hard to drain away.

Use PERFACTION few minutes a day it is enough to discover the advantages that join wellness and beauty.

This manual is suitable to quickly familiarize with your PERFACTION, permitting you to make a more safe and suitable use of it and to correctly effect the little maintenance operations that guaranty safety, reliability and length.

### **WARNING**

The performance of any operation with this appliance other than those recommended in this instruction manual, may present a risk of personal injury: possible personal injuries and material damages (appliance included) that led to, could be not ascribed to the manufacturer.

Power supply

- Before connecting the appliance to the electrical grid, always check if the mains voltage corresponds to the voltage indicated on the rating plate on the bottom of the appliance, and that has been equipped with safety devices in full working order.
- Before connecting/disconnecting the cable plug to an earthed power supply with correct voltage, make sure that the on/off switch is on "off".
- Ignition and turning off must be effected by the switch (for ex. do not unplug the power supply cable when the appliance is working)
- If the appliance stay unused for a long while, always disconnect it from the mains.

Warning!

Do not use the appliance with wet hands or if some water is near

### **Placing and Protection**

- Do not expose for a long-time the appliance to direct sunlight, to heat source, to strong vibrations, to very wet and dusty places, to the rain, because all them could be cause malfunctions, deteriorations and electric shocks also.
- Storing the appliance always turned towards to the up side, with the handle inserted into its lodging and with the tube re-winded into the proper place.
- Be careful to not damage your heartbeat detector. Always use it with extreme care and put it in a safe place.

### **Duty supervision**

PERFACTION is equipped with safety measures to safeguard the users. However every appliance may present

so-called “non-eliminable residual risks”, or rather non evident potential risks that can’t be totally reduced with planning and safety systems.

Before using the appliance is compulsory read these instructions and all the attached documentation, giving special attention to the prescriptions concerning the safety.

Always watch over that the appliance is not utilized by people haven’t got the requirements to do it, cause age or psychophysical capabilities (as specified as follow). Particularly every form and purpose of use of the appliance (also different from the intended use), must be forbidden to children.

When not in use and/or non supervised, the power supply cable must be unplugged from the mains, and the heartbeat detector must be disconnected also.

### Other use restrictions

The appliance use can be effected only from:

- Adults with all the common psychophysical capabilities.
- Under ages from 14 years old with all the common psychophysical capabilities, even only under parents or responsible people supervision but whatever after medical advice only.

The appliance use is explicitly prohibited to:

- Children under 14 years old
- People with reduced common psychophysical capabilities.

NOTE: some particular states of physical tiredness and/or debilitation, also if not caused by working activity or by clear and/or certified illness, must be considered and treated as absence of all the common psychophysical capabilities.

- People under influence of alcoholic drinks, tranquilizers or any kind of drugs.

NOTE: in these cases must be strictly prohibited approaching to the appliance also.

### Warning!

People use and/or give in an accommodation the above described appliance, explicitly acknowledge that none responsibility can be attributed to the manufacturer in case of accident occurred in presence (and not necessarily “in consequence”) of violation of one or more dispositions (none excluded) indicated on this document.

### ADVICES

#### Warning!

- Do not use this appliance unlike as recommended in this instructions manual, on the packing and/or in the relative illustrations, it could be cause of risks for the users’ safety.
- The manufacturer and the sales network decline their liability for any damage or accident derived from any use not in conformity with the instructions contained in this manual and/or in the warnings and/or on the illustrations.
- Always remember that the non respected instructions indicated on this manual violate the legislation in force and, in case of accident, it could have a criminal grave at your burden.
- Non respected instructions indicated on this manual cause the lapse of every kind of guarantee definite and/or implicit.
- All the documents attached to this manual are an integral part of it, and the reading of them is necessary and can’t be fail.

Warning! Never attempt to modify the appliance:

- The appliance manufacturer and the retailer are not responsible for any direct or indirect consequence for yours or other people’s alterations.

- Always remember that the appliance alteration is a violation of the legislation in force.
- Every appliance alteration cause the lapse of every kind of definite and/or implicit guarantee.

### **Cleaning and ordinary maintenance**

A regularly cleaning guarantee the working order, the good functional and the long life of the appliance.

Particularly:

- Make sure that the appliance filters, situated in the hollow under the handle, is clean from impurities; wash and dry it with care periodically.
- Avoid the sediment formation and that some dirt get into both the lodgings of the suction tube and the heartbeat detector, and into the handle suction holes.

No other maintenance operations are requested.

### **Warning!**

Before attempting any maintenance operation, always disconnect the appliance from the mains and disconnect the heartbeat detector.

To clean the body of the appliance, simply use a soft dry or a moist cloth.

Do not use: alcohol, benzene, acetone, solvents, detergents, abrasive products, chlorinated substances, steel-wools, steel or synthetics brushes, scrapers, etc.

### **Warning!**

Do not immerse totally or partially the appliance and its parts into water or other liquid fluids or utilize water in jets or any other origin and quantity.

### **Repair and special maintenance (Only for Italy)**

The special maintenance and the repairs are very delicate operations, good experience is required and must be effected from skilled workers and entrusted directly from the manufacturer.

In case of service necessity:

- The end user can contact directly and promptly to give indications on the number 199,151,136 authorized repair center to open the practice of care.
- The retailer can start the standard of care Joycare requesting authorization to return by e-mail to: assistenza@joycare.it
- Never use the appliance if it's in full working order.
- Do not open the appliance body.
- Never attempt to replace and/or repair any parts of the appliance.
- Never try to solder and/or stick damaged parts.

### **Manufacturer guarantee**

- This product is intended for household use only and not for professional use, it's intended directly to the customer, and is guaranteed 24 months from the purchase date.
- The guarantee assures the free replacement of parts damaged or not fair wear and tear caused by a checked defect of planning, supply, construction and/or assembling to attribute to the manufacturer.
- For repairs caused by a defect in workmanship, the charges for the labour are at manufacturer's care, within the limits to interventions of technical personnel entrusted by the manufacturer, in application of this guarantee.
- All the forwards for the dispatch and the restoration of the appliance to repair, are made free during the guarantee period and after. For more information call **199.151.136**

**Important:** in accordance with the Italian legislation in force, that acknowledges and brings about EU directives, the manufacturer is the one who affix one's marks or other distinguishing signs on the product or that presents himself as the manufacturer: therefore the manufacturer doesn't respond for products lacking in the above mentioned marks and distinguished signs and/or tempered and/or covered and/or disguised.

**Warning! The manufacturer guarantee does not include:**

- Appliance replacement after a not authorized repair.
- Any refund for working hours effected from personnel non authorized directly from the manufacturer to make some repair.
- Any damages caused by the use of the appliance into not suited rooms.
- Any phenomena not due to regular working of the appliance (electrical atmospheric discharges, tension leaps, fires, etc.)
- Any refund for non-earnings (direct or indirect), and/or for expenditures (direct or indirect), and/or penalties (direct or indirect), and/or damages (direct or indirect) in business or in image caused by the appliance stop.

**Warning!** The manufacturer's guarantee totally lapse in the following cases:

- Violation of one of the some dispositions, instructions, recommendations, warnings and/or of one of the some manufacturer's warnings and/or prohibitions. Particularly in case of violation of one of the some dispositions (of using and not only) written on this manual and/or on any attached document.
- Any improper or unreasonable use of the appliance and/or its parts, even if not foreseen and/or not expressly prohibited by the manufacturer.
- Any alteration, slight also, of the appliance and/or its parts.
- Any replacement of appliance parts and/or components not resulting from the manufacturer, independently from their technical characteristics.
- Any repair or technical intervention of personnel not authorized from the manufacturer.
- Any appliance damage given by heedlessness and/or negligence.
- If the appliance wasn't duly purchased and paid to the manufacturer or the vendor, in accordance with the terms of the agreement. Particularly, every late effected or not-payment due to lapse the guarantee validity of both the manufacturer and the vendor during all this period. These periods can't be recovered after debt payment. Whatever the manufacturer's guarantee lapses 24 (twenty-four) months after the delivery date.
- Without regular invoice or receipt or payment check showing the regular purchase.

**What for using PERFACTION**

Listen to your heart means take care of yourselves, respecting the organism and his needs with effective gestures, without looking for some "miracles".

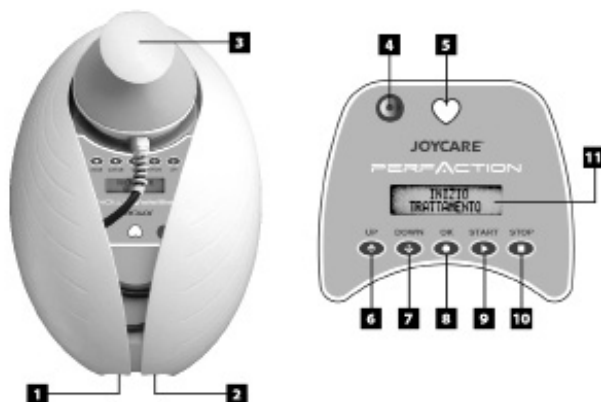
Swan listen to your heartbeat, acting in synchronicity with it, adapting his action to your body and offering great efficacy and comfort to you.

Use it regularly in all the following cases:

- Reducing the located adiposities and free yourselves from the cellulite.
- Smoothing down the peel skin and the stretch marks.
- Stimulating blood circulation and lymph draining.
- Reducing the thighs measurement, toning up gluteus and returning tonicity to the tissues.
- Feeling good with your body.

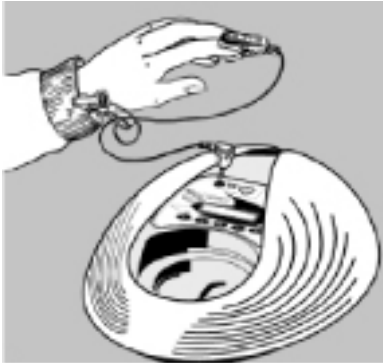
## Ignition and utilization

- 1 SWITCH
- 2 POWER SUPPLY
- 3 HANDLE
- 4 HEARTBEAT SENSOR SOCKET
- 5 HEARTBEAT LIGHT
- 6 UP•TIME+ KEY
- 7 DOWN•TIME- KEY
- 8 OK•CONFIRM KEY
- 9 START KEY
- 10 STOP/PAUSE KEY
- 11 DISPLAY

- Plug the power supply cable into its site 2 and connect it to the mains.
- Plug into the socket 4 the heartbeat detector connector and apply it on a finger of the hand that doesn't hold the handle, fixing it on the finger tip and on the nail.
- Unlock the handle 3 from its lodge and unwind the tube.
- Switch on the appliance 1
- Wait for confirmation of heartbeat detector, indicated by an **acoustic signal**
- Set on the display 11 the time for the required treatment, using up & down bottom 6 and 7
- Press ok 8 to confirm the time treatment setting.
- Sit down comfortable and push **start** 9 to beginning the treatment. The **heart-shaped** illuminated sign 5 blink
- When the treatment setting time is over, the appliance automatically stop, emitting an acoustic signal. Then pushing pause/stop 10 you can set a new treatment.
- By pushing pause/stop 10 during the treatment :  
One time: the treatment goes to pause mode;  
 By pushing start 10 the treatment restarts from that point;  
Two times: the treatment was definitely stopped (now it is possible to set a new treatment)

## Heartbeat detector

The **PERFACTION** secret is kept in your heart.

The heartbeat detector sends to the **PERFACTION** the signal that controls the working of it.

Thanks to this patented technology the **PERFACTION** working is always synchronized to your real personal condition, minute after minute.

So you can obtain greatest results, always in respect of your life cycles.

Apply the detector on a finger of the hand that doesn't use the handle, fixing on the finger tip and on the nail.

To obtain a good and more exact reading of the heartbeat, keep still hand during the treatment.

NOTE: During the treatment, the display could signal a check request of the heartbeat detector. It can depend on a reading defect due to the losing touch with the detector or to its uncorrected position.

Provide to correctly re-positioning it and go on with it.

In case of persistent heartbeat reading shortage, the appliance automatically stops the treatment until a good reading will be restored.

NOTE: If during the treatment you don't feel the suction effect from the handle on the skin (VACUUM effect), control that the closing filter do not leaking air.

### Suggestions and use cautions

You can use PERFACTION on all the external whole parts of your body (except for eyes, ears, glottis, nipples and genitals), with the following cautions:

- Execute the treatment on dry and clean skin to facilitate the handle sliding.
- During the treatment apply a definite but not excessive pressure on the handle: the use of PERFACTION never be cause of pain or persistent blushing or impress any marks on the skin. Contrary, interrupt the treatment and ask your doctor's advice.
- Avoid brusque or tiring movements, that could caused heartbeat sudden changes.
- If you feel an excessive vacuum effect, slightly lift one of the sides of the handle to nullify it.
- To obtain a good result, begin making the treatment at least three times a week, for three consecutive weeks.

After that, continue with at least one/two times a week as maintaining action.

The PERFACTION never been used :

- Into the oral hollow.
- On the mucus membranes

- On the atypical appearance moles, only after doctor's advice.
- On jewels or piercing present in the skin areas to treat.

The PERFACTION never been use peremptory in the following cases:

- Children under 14 years old (for minors over 14 years old ask your doctor's advice).
- Women in pregnancy.
- People suffer from epilepsy, mental or behaviour troubles.

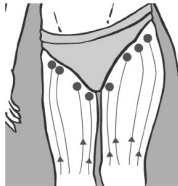
Furthermore PERFACTION never been used in the following cases: fever, infections, skin diseases, inflammations, open or recently healing wounds, blood coagulations troubles, pressure serious troubles, phlebitis, surgical operations to the vascular system, varices and/or pins and needles to the limbs. To clear up any doubt ask your doctor's advice.

In case of prolonged consumption of medicines as antibiotics, antihistamines, anti-inflammatory and anticoagulants before utilizing PERFACTION ask your doctor's advice.

## WOMAN

### Cellulite program

#### Thighs: 15 min.

Starting from the knee, move the handle very slowly from down to up linger over where the slight blemishes are more in evidence.

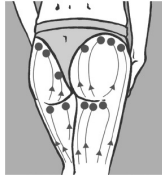
Insist for few seconds on the indicated points (●)

Repeat the treatment on the other leg.

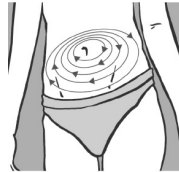
#### Belly: 15 min.

Massage very slowly around the navel, proceeding with a circular movement in clockwise direction.

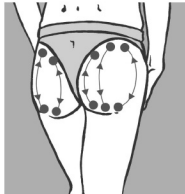
**Thighs / Gluteus: 15 min.**

Starting from the knee move the handle very slowly from down to up, linger over more evident slight blemishes. Insist for few seconds on everyone of the indicated points (●) Repeating the treatment on the other leg.

**Belly / Hips reduction**

**Duration: 15 min.**

Massage with a round and very wide movement to treat the more external parts, without fear to feel a more accentuated vacuum effect.

**Hips toning up**

**Duration: 15 min**

Start from the indicated points (●) in the basis of the gluteus.

Massage in definite way, following the muscular line of the gluteus and moving it with the handle, without fear to feel a more accentuated vacuum effect.

Insist for a few seconds on everyone of the indicated points (●), proceeding then to the next one

**Legs treatment**

**Duration: 15 min.**

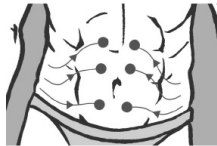
Stimulating action of the micro-circulation, very useful to maintain tonicity and resilience of the tissue.

Starting from the knee move the handle very slowly from down to up, insist for few seconds on the indicated

points (•).

Repeat the treatment on the other leg..

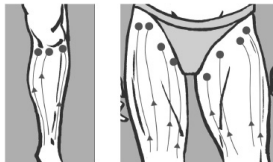
## MAN

### Abdominals / Hips remodelling action

**Duration: 15 min.**

Massage starting from the hips and move to the navel (as illustrated), insist for few seconds on the indicated points (•)

### Decongestant thigh / calf length: 15 min.

Starting from the calf, very slowly move the handpiece from the bottom up, pausing for several seconds on the discharge points indicated (•) Repeat the treatment on the other leg.

### Action Remodeling pectoral

**duration: 15 min.**

Massage movements semicircular (as shown), insisting on for a few seconds discharge points indicated (•)

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