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Operating instructions

BS37 is a multifunctional 6-in-1 infrared forehead and ear thermometer, which lets you check the exact body temperature of your baby or an ailing person with an easy and uncomplicated method.

The device also has a scan and stopwatch function. The device can moreover measure the room temperature and also display the date and time of day.


1. Scope of delivery

1 x BS 37 forehead and ear-thermometer
1 x battery CR2032
1 x operating instructions

2. Features

- infrared measurement of body temperature on the forehead and in the ear
- precise measured values in a few seconds
- alert at temperatures over 37,5 °C/99,5 °F
- stopwatch function
- storage of the last 9 measured values
- 12/24-hour display
- switchover from °C to °F

3. Control elements

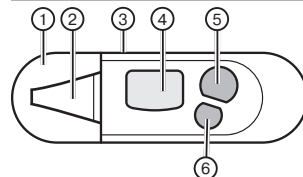1. Protective cap
2. Infrared sensor
3. Scan key for ear measurement
4. LCD display
5. Scan key for ear measurement
6. Multifunctional key

4. Safety notes

Any warranty claim is void in the case of damages caused by non-observance of these operating instructions. We are not liable for consequential damages or for damages to property or persons caused by

non-observance of the safety notes and improper operation of the thermometer.

- Please read the operating instructions before using the thermometer!
- Arbitrary modifications and/or alterations of the product are not permitted for safety and approval reasons (CE).
- BS 37 may only be operated with a battery (1 x button cell CR2032). Never try to operate the device with a different power supply.
- The device and the batteries do not belong in the hands of children – they are not a toy. Batteries can be lethal, if they are swallowed. Immediately obtain medical help, if a battery has been swallowed.
- Immediately switch the device off, if it:
 - has visible damages,
 - does not work anymore.
- Only qualified specialist staff may carry out repairs and adjustment work.
- Avoid the following influences to the device:
 - strong mechanical stress,
 - extreme temperatures,
 - strong vibrations,
 - high humidity.
- Do not leave the packing material lying around. Plastic films, plastic bags, cardboard boxes etc. could become a danger-ous toy for children.

5. Useful notes

- The body temperature of the human body varies during the day. The lowest temperature is in the early morning and the highest in the evening.
- You should take your temperature and that of your family members when you/they are not sick and note the measured value. You can then immediately recognise when measuring in the case of sickness, whether the sick person is running a temperature.

6. Functional overview


Time/date: The time/date display helps you read the accomplished measurement in combination with the date and time memory function. The time/date display changes automatically every 5 seconds on the display in

the OFF mode. (See section 7.1. *Time/date setting* to carry out the required settings.)


Time:

Date:

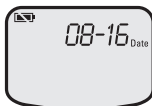Room temperature: The current room temperature is always displayed in the OFF mode. This is recognisable by the symbol to the right of the temperature display (see section 8.3. *Measurement of room temperature*).

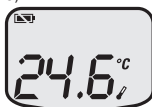Ear or forehead measurement: The forehead measurement is meant to determine illnesses in combination with a temperature – it never replaces consulting a physician! (See sections 7.2. *Use of the ear-thermometer* and 7.3. *Use of the forehead-thermometer*, to carry out the appropriate settings.)

Scan function: You can measure surface temperatures in the scan mode, e. g. the temperature of milk. (See section 8.1. *Scanning*, to carry out the appropriate settings.)

Stopwatch function: The stopwatch function is particularly suitable for measuring the heart beat (see section 8.2. *Stopwatch*).

Temperature alarm: A long beep, followed by three short beeps acoustically confirms that the BS 37 has registered a body temperature of more than 37.5 °C/99.5 °F. This indicates a high temperature.

Memory function: BS 37 has 9 memory locations, so that the last 9 measurements can be stored and called up again. You are in the memory mode, if the appropriate symbol is visible on the display below the battery symbol. (See section 8.4. *Measured value storage*, to carry out the appropriate settings.)

Switching function °C/°F: You have the option of measuring the temperature in °Celsius or in °Fahrenheit (see section 8.5. *Switching from °C to °F*).

7. Start-up

7.1. Time/date setting

1. Press the multi-function key (6) until a double beep is heard.
2. Then press and hold the multi-function key (6) until the "SET" symbol flashes on the top right of the display.
3. As soon as the "SET" symbol flashes on the display, briefly press the scan key (ear or forehead measurement), so that the number "24" flashes at the top left of the display. The time is then displayed in the 24-hour mode. If you decide for the 12-hour mode, press the scan key (ear or forehead measurement) again, so that "PM" appears (after 12 o'clock) or no symbol appears (before 12 o'clock).
4. You can now enter the current hour, by briefly pressing the multi-function key (6) once after the 12/24-hour input, so that the hour display flashes. Now enter the appropriate hour with the scan key (ear or forehead measurement).
5. To enter the minutes and date, press on the multi-function key (6) again, so that the appropriate minute and/or date display flashes. Enter the appropriate data (minutes, year, month, day) by means of the scan key (ear or forehead measurement).


Note! The device automatically switches to normal operation, if no changes are made in the setting mode within five seconds.

7.2. Use of the ear-thermometer

1. Ensure prior to use of the thermometer that the infrared sensor (2) is clean and dust free. To activate BS 37, press the multifunction key (6) once so that the display is shown as in the diagram for forehead/ear measurement.

4. As you are pulling the ear, carefully place the measuring tip into the ear as far as possible to close off the ear channel.

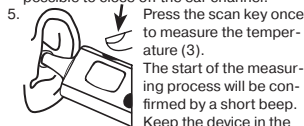ear until the thermometer beeps twice to signal the completion of the measuring process. Now you will be able to read the measured temperature on the LCD display.

When measuring temperatures of children under the age of one, pull the ear straight to the back.

When measuring temperatures of children over the age of one, pull the ear up diagonally.

- The measurement is normally accomplished within 3 to 8 seconds, depending on how much time the BS 37 requires to determine the body temperature.
- In some instance the measured temperature in the right ear, may differ slightly from the measured temperature in the left ear. For this always carry out repeat measurements in the same ear.
- Laying on the ear that is then taken for measurement may increase the measurement results. In this case, we recommend waiting a few minutes before measuring the temperature.

Important!

Slight measurement fluctuations ($\pm 0.2-0.3\text{ }^{\circ}\text{C}/\pm 0.4-0.5\text{ }^{\circ}\text{F}$) can occur, as is the case for all other thermometers, if several sequential measurements are made. In order to reduce measuring errors, always allow 30 seconds to elapse between individual measurements.

7.3. Use of the forehead-thermometer

1. Ensure prior to use of the thermometer that the infrared sensor (2) is clean and dust free. To activate BS 37, press the multifunction key (6) once so that the display is shown as in the diagram for forehead/ear measurement.

Hold the device to the temple and press the scan key once for the forehead measurement (5). The start of the measuring process will be confirmed by a short beep.

3. Hold the device to the forehead until a two-fold acknowledgement sound signals the end of the measuring procedure. You can now read the measured temperature on the LCD display.

Note!

- The symbol to the right of the temperature display flashes continuously during the measuring procedure. The symbol °C is constantly visible in the display after completion of the measurement.
- The measurement is normally accomplished within 3 to 8 seconds, depending

on how much time the BS 37 requires to determine the body temperature.

- Both the left or right side of the forehead can be used for the forehead measurement.

Important!

Slight measurement fluctuations ($\pm 0.2-0.3\text{ }^{\circ}\text{C}/\pm 0.4-0.5\text{ }^{\circ}\text{F}$) can occur, as is the case for all other thermometers, if several sequential measurements are made. In order to reduce measuring errors, always allow 30 seconds to elapse between individual measurements.

8. Functions

8.1. Scanning

1. The device must be turned on to activate the scan function.
2. Press and hold the multi-function key (6) when the device is turned on. Simultaneously press the scan key (ear or forehead measurement) briefly once, so that the symbol °C is also displayed at the bottom right of the LCD display.
3. To measure the temperature press and hold the scan key (ear or forehead measurement) for the duration of the examination. You can now e. g. check the temperature on a warmed up bottle of milk or a glass with baby food. The measured temperature is then displayed on the LCD display.

8.2. Stopwatch

1. The device must be turned on to activate the stopwatch function.
2. Press and hold the multi-function key (6) when the device is switched on, simultaneously press the scan key (ear or forehead measurement) twice successively, so that the symbol °C is displayed at the top left of the LCD display.
3. The stopwatch function starts, if you now press the scan key (ear or forehead measurement) in this mode. Renewed pressing of the key (ear or forehead measurement) interrupts the stopwatch function. Actuating the key (6) resets the stopwatch to zero. The stopwatch works for a period of 0.01 seconds to 3 minutes.

Example heartbeat measurement: Start the stopwatch and measure e. g. someone's pulse for 15 seconds. If you measured 18 heartbeats in 15 seconds, multiply the value by 4. The result is therefore 72 heartbeats per minute.

Note! The device automatically switches to off approx. 1 minute after the last measurement.

8.3. Measurement of the room temperature

The device automatically measures the room temperature as soon as the device is in the OFF mode. The room temperature is displayed by the symbol °C under the temperature display. The room temperature display updates itself automatically; this can take up to 15 minutes.

Note! The device should not be exposed to direct sunlight to determine the room temperature.

8.4. Measured value storage

The device automatically stores the last measured 9 temperature values. The value measured last is thereby always allocated to position 1. In order to call up the memory, press the multi-function key (6) when the device is switched on, so that the symbol °C is displayed top left on the LCD display. During the measuring, the symbol for ear measurement °C or the symbol for forehead measurement °C will be shown on the right side of the display. The memory position number always appears first followed automatically by the measured temperature. Renewed pressing of the multi-function key (6) displays the next memory position number.

Note! The temperature values are only stored in the memory, if the device is able to automatically switch to the OFF mode after the measurement.

8.5. Switchover from $^{\circ}\text{C}$ to $^{\circ}\text{F}$

To switch from $^{\circ}\text{C}$ to $^{\circ}\text{F}$ hold down the scan key for the ear measurement in the Off Mode (room temperature display) and press the

multifunction key (6) for approx. 3 seconds, so that the display will show the 'F' symbol. Follow the same procedure to switch the temperature from 'F' to 'C'.

9. Replacing the batteries

Replacing the batteries is necessary as soon as the symbol  appears on the display. For this purpose open the compartment on the back of the device and remove the batteries from the holder by means of a lever-age tool (e. g. screwdriver). Now insert a new battery (Type CR2032), while observing the correct polarity, into the device and close the battery compartment cover.

10. LCD display – Error display

Error display	Problem	Solution
	Error display 5-9 – the device is not working correctly.	The battery is drained. Insert a new battery into the device.
	The measuring procedure was started too soon.	Only proceed with the measurement when all symbols can be lastingly seen on the display.
	The device is subject to very fast changes of temperature.	Put the BS 37 into a room with a room temperature of 10 °C to 40 °C (50 °F to 104 °F) for approx. 30 minutes.
	The ambient temperature is not in the permissible range of the device (10 °C to 40 °C/50 °F to 104 °F).	

Error display	Problem	Solution
	1. The measured temperature exceeds 42.2 °C (108 °F) in the forehead/ear modus. 2. The measured temperature exceeds 80 °C (176 °F) in the scan modus.	Please note that the to be measured temperature must be within the limit values of the device.
	1. The measured temperature is less than 34 °C (93.2 °F) in the forehead/ear modus. 2. The measured temperature is less than -22 °C (-7.6 °F) in the scan modus.	
	The device cannot be activated.	Please replace the battery..

11. Notes on Disposal

 Old units, marked with the symbol as illustrated, may not be disposed of in the household rubbish.

    Used rechargeable and non-rechargeable batteries which are marked with one of the symbols illustrated may not be disposed of in the household rubbish.

You must take them to a collection point for old units, old batteries or special waste (enquire at your local authority) or the dealer from whom you bought them. These agencies will ensure environmentally friendly disposal.

12. Care and warranty

Before cleaning the unit, disconnect it if necessary from other components; do not use aggressive cleaning agents. To guarantee an accurate measurement, use alcohol to clean the sensor prior to each measurement. The unit has been carefully checked for defects. If nevertheless you do have cause for

complaint, please send us the unit with your proof of purchase. We offer a 2 year warranty from date of purchase.

We are not liable for damage arising from incorrect handling, improper use or wear and tear.

We reserve the right to make technical modifications.

13. Technical data

Ambient temperature range:
10 °C-40 °C (50 °F-104 °F)

Forehead/ear measurement:
34 °C-42.2 °C (93.2 °F-108 °F)
Forehead/ear measurement tolerance:
± 0.2 °C (0.4 °F) at 35.5 °C-42 °C (95.9 °F-107.6 °F) for ear measurement and at 36 °C-39 °C (96.8 °F-102 °F) for forehead measurement; ± 0.3 °C (0.5 °F) for all other measurement ranges

Scan mode:
-22 °C-80 °C (-7.6 °F-176 °F)
Scan mode tolerance:
± 0.3 °C (0.5 °F) at 22 °C-42.2 °C (71.6 °F-108 °F); ± 4 % or ± 2 °C (4 °F) for higher temperatures

Display:
Degrees (°C) or Fahrenheit (°F)

Memory:
the last 9 measurement results are automatically stored

Battery type:
CR2032

Battery life:
more than 1000 measurements

Automatic deactivation:
The device switches itself off automatically 1 minute after the last measurement.

Dimensions:
123 mm x 37 mm x 30 mm (l x w x d)

Weight including batteries:
62 g

This device complies with the regulations of the EEC Directive 93/42/EEC.

 0197

You can find up-to-date product information on our website <http://www.hartig-helling.de>

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