innofit

easy tech, easy life

DOLCE TEPORE

INN-051
MANUALE D'USO

USER MANUAL

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Congratulations on having purchased Innofit DOLCE TEPORE model INN-051, ready-to-use hot water bag. Before starting using the appliance, please read carefully all the instructions below. It is advisable to save this instruction manual for a future use.

IMPORTANT SAFETY INSTRUCTIONS

Never use the product when it's connected to the electric power (heating phase). When it's heated, please disconnect it from the power supply, it's ready to use.

Before starting using the appliance, please check the product is not damaged. In case of doubts, please address to a specialised centre.

In case of damage, functioning problems, or falling into water, do not use the product, please contact the after sales service. Always pulling from the plug and not from the cable. Do not pierce the fabric of the product and before use, always check the product has no leakage. Do not swallow the liquid of the bag. Never try to change the water inside the bag. Before connecting the appliance, check the voltage shown on the rating label corresponds to the local mains voltage. Do not place or store product where it can fall or be pulled into a tub or sink. Do not place it or drop it into water or other liquid. In case the product falls into water **NEVER TRY TO REACH IT!**

Always store the product away from heating sources, sunlight, dampness and sharp

objects. Never use the product while having a bath or a shower. Never touch the product or its power cord with wet hands. Never place the product on a wet or hot surface, this might damage the unit. In order to protect children from electrical shocks or other accidents, never leave the appliance unattended when it is plugged in and always store it out of their reach. Disabled and people insensitive to heat must not use this product. Use this product only for its intended purpose, as described in this instruction manual. Any other use is improper and therefore hazardous. Use original accessories only. The product is for an indoor use only. Servicing of this product should be done by qualified centres only and replacing components must be original.

HOW THE PRODUCT WORKS

Take out the water bag from the pouch and place the bag with the plug upwards on a flat and stable surface. Plug in the power cord in the plastic part of the bag and connect the plug to socket. DOLCE TEPORE heats up rapidly (10-15 minutes) and safely. At the end of the heating-phase, the temperature of its surface reaches approximately 50°C and the LED indicators turns off. The heating time is about 15 minutes. Disconnect the power cord from the socket and the product. The hot water bag is ready to use. INN-051 remains hot for 2 to 6 hours, depending on the outside temperature.

When the bag is cool, it's possible to heat it up again.


After the use, always put the bag inside its pouch and store it in a clean and dry place, keep it away from heating sources and sun rays. Humidity might damage the product.

In case of skin eruptions, wounds, inflammations, blood circulation ailments address to your physician before treating this part with the product.

This product is not a medical device and it's not suitable for professional medical treatments. This appliance is for home use only.

In case of rheumatic pains, dysmenorrhoea, cephalaea, muscle/bone pains or any other diseases, please consult your doctor before use.

RATING LABEL:

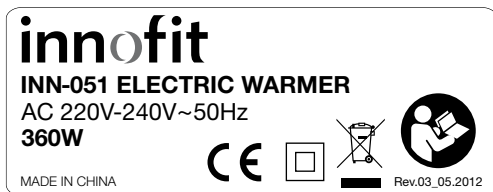■ This appliance is placed on the market after 13th August 2005 and must be separately collected.

Separate Collection- WEEE (Waste of Electric and Electronic equipment)

At the end of its use don't dispose the appliance, including removable parts and accessories, together with the other urban waste, but conformably to 2002/96/EC. Since it has to be treated apart from home waste, you must take it to a differential collection centre specific for electric and electronic equipment (EEE). Otherwise you can give it back to the retailer as you buy an equivalent apparatus.

There will be severe sanctions in case of transgression. The correct disposal helps preventing potentially negative consequences on the environment and human health.

Read carefully the instruction manual.

This product complies with all the applicable EU directives.

NOTE: Specifications and designs are based on the latest information available at the time of printing and subject to change without notice.

innofit

easy tech, easy life

REV.03_05.2012

Innofit S.p.a.

Via Merloni, 9 • 60131 Ancona • ITALY • TEL +39 071 2867616 • FAX +39 071 888165
www.innofit.it