AUDIOLINE

Bedienungsanleitung
Operating Instructions
Mode d'emploi
Istruzioni d'uso

Watch&Care FS-700

D

GB

FR

IT

⚠ WICHTIG

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IMPORTANT

- Measure the temple area or behind the ear lobe, carefully following the instructions.

Temple area

or

Behind the ear lobe

Incorrect measurement

If there is sweat on the forehead, measurement from behind the ear lobe is recommended.

2 ~ 3 cm

2-3 cm away from the forehead.

Wipe off sweat

Push back hair

Wipe off any sweat and push back any hair from the forehead before taking a measurement.

Do not take any temperature measurements when near to a stove or air conditioning unit. Do not use outdoors.

Wait 30 min.

If there is a significant change in ambient temperature, allow both the thermometer and the patient to adjust for at least 30 minutes. Wait 30 min. before taking the patient's temperature.

1	General Information	22
	Brief product introduction	22
	Basic operation principles	22
	Product features	23
	User precautions	23
	About body temperature	23
2	Operating Elements	27
3	Getting started	28
	Batteries	28
	Switching on and setting the thermometer	28
	Setting the measurement unit	29
	Setting the measurement sound	29
4	How to use FS-700	30
	Measuring body temperature	30
	Displaying stored measurements	32
	Measuring object temperature	32
5	Appendix	34
	Error messages	34
	Maintenance and Storage	35
	Product Specifications	36
	Disposal	38
	Guarantee	38

General Information

1 General Information

Please read these operating instructions carefully. If you should have any questions, please contact our Hotline. Please refer to our website **www.audioline.de** for the telephone number. Under no circumstances open the unit yourself and do not try to carry out your own repairs. For any claims under the warranty, please refer to your dealer.

Brief product introduction

This manual is the user manual for the IR Thermometer (Model: FS-700). This product is a Non-Contact IR Thermometer.

- The FS-700 takes contactless temperature measurements, ensuring that the body temperature can be measured safely.
- The contactless design prevents the possibility of infection transmission associated with normal thermometers.
- The FS-700 can measure the body temperature safely and without contact.
- The FS-700 measures the temperature in less than 2 seconds.
- Please read the manual carefully.

Basic operation principles

- When your body is hot, you touch your forehead because body's temperature is best reflected there.
- The thermometer measures infrared energy radiating from the forehead. Once the thermometer has analysed the amount of infrared energy via a sensor, it converts it into a temperature value.
- The FS-700 is a contactless thermometer which allows the user to measure the body temperature about 2 to 3 cm away from the forehead, thereby decreasing the possibility of infection or discomfort resulting from contact.

Product features

- LCD light for ease of use and digital LCD display.
- Multi-measurement function which measures both body and object temperatures.
- Storage of up to 32 results.
- Contactless measurement for increased safety.
- Power-saving: Automatic shut-off after one minute.
- Powered by two AAA Type Alkaline batteries.

User precautions

- Use this product from 16 °C to 40 °C.
- Do not expose the product to temperatures below -20 °C or above 50 °C or above humidity levels of 95% ($\geq 95\%RH$).
- Measurements taken using this product do not represent a doctor's consultation or examination.
- Examine the product periodically for any damage/faults.
- Read the manual carefully before using this product for the first time.
- In the case of any product damage or malfunction, do not use it and request its repair through the authorized dealer.
- If the thermometer is not used for a long time, please take out the batteries and store it in a safe place. The used batteries must be collected separately.

About body temperature

- What is body temperature?

The body temperature is the human body's internal heat level. The 'normal' body temperature is that of a healthy person and can vary slightly within certain limits. The body temperature changes according to the environment, time, and activities that a person might be doing. For example, a person's body temperature can often be 0.5 °C higher at night time than early in the morning.

General Information

- Body temperature control

The central nerve in the brain's hypothalamus controls the body temperature, maintaining it at a certain level. When the body temperature goes up, it expands the blood vessels inside the skin and releases sweat from the sweat glands, resulting in significant heat loss – a hypothalamic reaction to suppress any further heat increase. However, when the body temperature goes down, it contracts the blood vessels and prevents sweat secretion from reacting to the stadium incrementi by reducing heat. The human body can maintain a relatively regular body temperature because heat production and loss always equates.

- Body Temperature Measurement Ranges


Normal ranges by area:

Axillary: 34.7 ~ 37.3 °C

Oral: 35.5 ~ 37.5 °C

Rectal: 36.6 ~ 38.0 °C

FS-700 (Forehead) : 36.0 ~ 37.9 °C

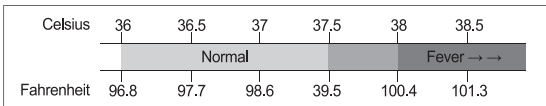- Why do fevers occur?

Normally, when bacteria, viruses or other pathogens attack the body, fever occurs. Other reasons for fever can be tissue damage or blood diseases. When burns and heart attacks occur, the blood vessels also contract and prevent sweat secretion from reacting to the stadium incrementi by reducing heat. The human body temperature can maintain a relatively regular body temperature because heat production and loss always equates.

- Infant body temperature

Infants cannot sufficiently control their body temperature, the surrounding environment will affect their body temperature. If the indoor temperature is high, the body temperature increases, and vice-versa. It is therefore important to check an infant's normal body temperature on a regular basis.

- Normal temperatures can be checked by calculating the average value for several days at regular hours in the morning and evening. An infant's body temperature can vary by about 0.5 °C according to the surrounding temperature. For example, when a temperature reading for 8 A.M. is 36.7 °C and 8 P.M. is 37.3 °C the normal temperature range is 36.7 °C to 37.3 °C and the normal temperature is 37 °C. Fever occurs when the body temperature rises 1 °C compared to the normal temperature.
- It is recommended that the indoor ambient temperature be set to 24 °C to 26°C, with humidity of around 40% to 60%. In winter, thermometers and hygrometers should be placed at the height of an infant's living area.
- Infants should not be covered with swaddles and laid on heated floors during summer. If miliaria occurs, or if the body temperature increases too much it can lead to dehydration fever. In the worst cases, this can lead to coma or cramps.
- Bathing environments for infants must be maintained at 24 °C to 27 °C. Five minutes of bathing with water at 38 °C to 40 °C is recommended. Infants must be dried well, and any water loss can be compensated through barley tea or fruit juice. During winter, the recommended bathing time for infants is between 11 A.M. and 2 P.M. Drying in cold air can lead to colds or cause cracks in the cheeks or hands.

General Information

- Breast-fed infants should drinking breast milk at its regular temperature, i.e. body temperature. The recommended temperature to serve powdered milks at is 38 °C.
- Using FS-700
 - Measuring Body Temperature
Switch the FS-700 to body temperature measuring mode, and take the body temperature from the temple area with a 2 to 3 cm distance.
 - Measuring Bathing Water Temperature
Switch the FS-700 to the daily wide temperature measuring mode and take the water temperature by placing the probe 2 to 3 cm from the water. Be careful not to get water on the FS-700.
 - Measuring a Nursing Bottle
Switch the FS-700 to daily wide temperature measuring mode and take the temperature by placing the probe 2 to 3 cm from the nursing bottle.

2 Operating Elements

1 LCD display
2 Start button
3 Sound button

4 MODE button
5 MEM button
6 Strap position

7 Probe
8 Battery cover

Buttons

Function

START

Switching the device on.
Starting temperature measurement.

(🔊)

Setting the acoustic signal (On/Off).

MODE

Setting the measurement method.

MEM

Displaying stored measurements (max 32).

MODE + **MEM**

Setting the measurement units (°C/°F).

Display

Meaning

Body temperature mode.

Object temperature mode.

Memory mode.

Change the batteries immediately.

Getting started

3 Getting started

Batteries

Before first using this product, insert the batteries supplied. The device requires two AAA (LR03) batteries.

- Open the battery compartment.
Slide down the battery compartment lid on the back.
- Insert the batteries.
Make sure that the batteries are inserted correctly.
- Close battery compartment.

Switching on and setting the thermometer

- Press the Start button **START** or Mode button **MODE** briefly.
The device switches on.
After two seconds, the display shows that the device is ready for use.
The device always starts in body temperature mode

Setting the measurement unit

- Press the Mode button **MODE** and MEM button **MEM** simultaneously for at least 3 to 4 seconds.
The device switches to the other measurement unit.
Possible measurement units: Celsius °C and Fahrenheit °F.

Setting the measurement sound

- Press the Sound button **🔊** briefly.
Each time the button is pressed, the sound switches on or off.
If the display shows **On** after pressing the button, the device emits sound.
If the display shows **OFF** after pressing the button, the device does not emit sound.
This setting is not stored. The next time the device is switched on, the sound is reset to **On**.


How to use

4 How to use FS-700

Measuring body temperature

- Press the start button **START** or Mode button **MODE** briefly.
The device switches on. After two seconds, the display shows that the device is ready for use. The device always starts in body temperature mode.
- Hold the thermometer 2 to 3 cm from the measuring point.
The ideal measuring point is the temple area.
- Press the Start button **START** briefly.
Maintain the distance until you hear the measurement completion beep. If the thermometer is in silent mode, maintain the distance for approximately 2 seconds.

Normal:

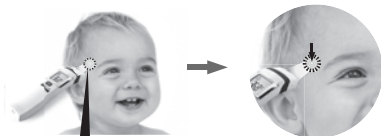Recommended measurement point is the temple area.

Sweating:

Behind the ear lobe.

Note

We recommend changing the measurement point to behind the ear lobe from the forehead and temple because this area is less affected by sweating than the forehead.

- Temperature readings may vary according to location. Therefore we recommend measuring the temperature more than two times on the forehead.

Recognize the highest temperature as the correct temperature.

- Check the temperature on the LCD screen. This will be automatically stored.
- The internal LED will be turned as the result shows up.
LED is automatically turned off after 2 to 3 seconds.

LED 2 colours

Green: measurement range 22 °C to 37.9 °C.

Red:   the temperature of the measuring point is above/below the measurement range.
 the device temperature is outside the operating temperature range 16 °C to 40 °C.

- When the measured temperature is more than 38 °C the fever guidance will help you to better explain the meaning of the infant's temperature.
The screen will be green if no fever.
- Read off the temperature on the display.
The measurement is automatically stored.

How to use

Displaying stored measurements

- Press the start button **START** or Mode button **MODE** briefly.
The device switches on. After two seconds, the display shows that the device is ready for use.
- The device automatically stores the last 32 measurements.
- Press the MEM button **MEM** briefly.
The most recent measurement is displayed.
The diagram shows two LCD displays. The first display shows a battery icon in the top right, a small 'r' on the left, and a vertical bar in the center. An arrow points to the second display, which shows the same battery icon, a small 'r' on the left, and the number '36.8' followed by a degree Celsius symbol on the right.
- Each additional press of the MEM button **MEM** calls up the next most recent measurement. The number on the display indicates how many measurements ago this reading was taken.
The thermometer then changes to body temperature mode after 30 seconds.
- The device only stores measurements when in body temperature mode.

Measuring object temperature

- Press the start button **START** or Mode button **MODE** briefly.
The device switches on. After two seconds, the display shows that the device is ready for use.
The device always starts in body temperature mode .
- Press the Mode button **MODE** briefly. The device switches to the object temperature measuring mode. .
- Check the measurement ready status on the LCD.

How to use

- Hold the thermometer 2 to 3 cm from the desired measuring point.
- Press the Start button **START** briefly.
Keep the distance until you hear the measurement completion beep.
If the thermometer is in silent mode, maintain the distance for approximately 2 seconds.

- Read the temperature from the display.
Measurements in object temperature mode are not stored.

Appendix

5 Appendix

Error messages

Q: The device temperature is outside the operating temperature range.

A: Operating temperature approx. 16 °C to 40 °C.

Q: The temperature of the measurement point is above / below the measurement range.

A: Measurement range for body temperature:
22 °C to 42.5 °C.

Measurement range for object temperature:
10 °C to 80 °C.

Q: The batteries are low.

A: Change batteries.

Q: Display appears briefly after changing the batteries.

Note

If you are exposed to temperatures of between 15 °C to 20 °C, the forehead temperature can decrease due to the influence of the ambient temperature.

Maintenance and Storage

- The sensor part of the thermometer is very sensitive and must be kept clean and undamaged.
- Do not put any foreign materials inside of the probe.
- To clean the sensor, use a cotton swab dipped in alcohol (ethanol) to carefully wipe the surface clean.
- The thermometer should be stored in a clean indoor environment, free from dust and pollutants and away from direct sunlight.
- The thermometer should be stored at room temperature. The suitable room temperature range is between 16 °C to 40 °C.
- If the thermometer is not stored at room temperature, leave the thermometer at room temperature for 30 minutes before using it.
- If the product needs be repaired, please contact the local authorized dealer.

Do not subject the thermometer to strong shocks, such as dropping it on the floor.

Do not submerge the thermometer or any of the components in water or other solutions.

Store the thermometer and the components in a clean, safe location.

Do not modify this thermometer without manufacturer authorization.

Do not disassemble or attempt to repair the unit or components.

Appendix

Product Specifications

Measurement method	Infrared (IR) measurement, Non-Contact measurement
Basic functions	Switches for multi-measurements (body and object temperatures) Switches for Celsius/Fahrenheit, Memory function, Sound controls (On,Off)
Operating conditions	16 °C to 40 °C, Humidity: below 95% Atmospheric Pressure (hPa) 700 ~ 1060
Transport and storage condition	-20 °C to 50 °C, Humidity: below 95% Atmospheric Pressure (hPa) 700 ~ 1060
Measuring distance	2 to 3 cm from the measuring point
Measuring range	Body temperature: 22 °C to 42.5 °C Accuracy: ± 0.2 °C (36 °C ~ 39 °C) Other Accuracy: ± 0.3 °C Object temperature: 10 °C to 80 °C Accuracy : ± 2 °C
Measuring time	Less than 2 seconds, continuous measurement
Resolution	0.1 °C
Mode selection	Body temperature, object temperature
Memory function	Automatic memory function Confirm the result by using the memory button. Can store up to 32 results
Sound	Button sound, measuring sound, result sound. The sound can be turned on and off by using sound button
Screen type	Liquid crystal display (mono LCD)

Display	Body temperature mode, object temperature mode, memory location, low Battery, °C/°F, Hi / Lo signs for the temperature outside the measurable ranges
Button	Start, Sound, Mode, Memory
Backlight	Turned on for 3 seconds after measurement The backlight turned off automatically
Power saving	Turns off automatically after one minute
Components	FS-700 1 set, Manual, Battery (AAA type) 2ea
Power	= 3V (AAA 1.5V Battery 2 EA)
Measurable numbers	More than 5,000 times
Size	152 mm (W) x 36 mm (D) x 37 mm (H)
Weight	68g (without batteries)

CE₀₁₂₀ The thermometer is in compliance with the following international regulatory and safety standards:
EN ISO 13485: 2012, ISO 9001: 2008, ISO 14971: 2012,
IEC 60601-1: 2005, EN12470-5: 2003 ASTM E 1965-98: 2003,
EN 60601-1-2: 2007, IEC 60601-1-6: 2010,
IEC 60601-1-11: 2010, IEC 62366: 2007, IEC 60601-1-6: 2010,
EN ISO 15223-1: 2012, EN 1041:2013, EN 62304: 2006

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Appendix

Disposal

In order to dispose of your device, take it to a collection point provided by your local public waste authorities (e.g. recycling centre). According to laws on the disposal of **electronic and electrical devices**, owners are obliged to dispose of old electronic and electrical devices in a separate waste container.

The adjacent symbol indicates that the device must not be disposed of in normal domestic waste!

Batteries represent a hazard to both health and the environment! Never open, damage or swallow batteries, nor allow them to get into the environment.

They may contain toxic and ecologically hazardous heavy metals. You are legally obliged to dispose of batteries at the point of sale or at appropriate collection points provided by your local public waste authorities. Disposal is free of charge. The adjacent symbols indicate that batteries must not be disposed of in normal domestic waste and but at authorised collection points.

Packaging materials must be disposed of according to local regulations.

Guarantee

AUDIOLINE equipment is produced and tested according to the latest production methods. The implementation of carefully chosen materials and highly developed technologies ensure trouble-free functioning and a long service life. The terms of guarantee do not apply to the batteries or power packs used in the products. The period of guarantee is 24 months from the date of purchase.

All deficiencies related to material or manufacturing errors within the period of guarantee will be redressed free of charge. Rights to claims under the

terms of guarantee are annulled following tampering by the purchaser or third parties. Damage caused as the result of improper handling or operation, normal wear and tear, incorrect positioning or storage, improper connection or installation or Acts of God and other external influences are excluded from the terms of guarantee. In the case of complaints, we reserve the right to repair defect parts, replace them or replace the entire device. Replaced parts or devices become our property. Rights to compensation in the case of damage are excluded where there is no evidence of intent or gross negligence by the manufacturer.

If your device shows signs of defect during the period of guarantee, please return it to the sales outlet in which you purchased the AUDIOLINE device together with the purchase receipt. All claims under the terms of guarantee in accordance with this agreement can only be asserted at the sales outlet. No claims under the terms of guarantee can be asserted after a period of two years from the date of purchase and hand-over of the product.

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Produttore: HuBDIC CO., LTD
Autore: Distributore: Inter KOTRA GmbH
AUDIOLINE GmbH