me™  
*smooth*

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## quick reference guide

### thank you!

Congratulations on purchasing your mē system! mē is an innovative device that can be used in the comfort and privacy of your own home to remove unwanted hair.

mē is the only at-home hair removal device to use patented elōs technology that is clinically proven to be safe and effective for all skin types and tones. Gentle, easy to use and just minutes a treatment, now you can experience the freedom of silky, smooth skin without the hassles of shaving and waxing.

Before using mē, visit our web site at [www.mepower.com](http://www.mepower.com) and make sure to watch our demonstration video.

### the most advanced hair removal technology

From Syneron Beauty and the scientists that invented IPL comes the new breakthrough patented hair removal technology called elōs. This unique combination of Intense Pulsed Light (IPL) and Radio Frequency (RF) addresses the limitations of traditional hair removal technologies and is clinically proven to safely deliver unprecedented speed and effectiveness for all skin tones.

### fast, easy and gentle

- **Fastest treatments.** Treat arms, legs, underarms and bikini area in just minutes.
- **Intuitive application.** Simply glide over your skin in a continuous motion similar to shaving. No tedious stamping.
- **Comfortable.** Just a gentle warming sensation – no nicks, cuts, stinging or snapping.

# quick reference guide

## how to use




### Step 1: preparing for mē

- Watch our demonstrational video at [www.mepower.com](http://www.mepower.com)
- Remove hair in the area you intend to treat.

### Step 2: sensitivity test

Before treating any new area with **mē**, test a small patch of skin near the area you intend to treat.

### Step 3: treating with mē

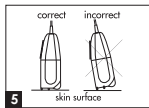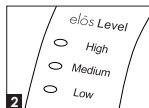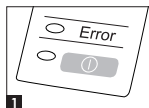1. Press the power button on the base unit to turn the system on.
  2. Press the appropriate elōs level button (Low, Medium or High) on the base unit to select the desired elōs energy. It is recommended to start with Low level, and if you feel comfortable, raise the elōs level.
  3. Press the elōs activation button on the applicator.
  4. Place the applicator on the area you want to treat, ensuring full contact of the two chrome (silver color) RF delivery bars with the skin. Do not press, just gently touch the skin. System will start emitting elōs pulses.
  5. Glide the applicator slowly and continuously on the treatment area without stopping in one place.
- Do not** return to the same spot for at least 10 seconds.

## is mē for you?

mē is intended for body use on ALL skin tones (skin color you are born with).

## quick reference guide

### body areas

mē is intended for use by women and men on body hair below the neck.

### get your professionally smooth skin

For best results, use mē on the same area once a week for 7 weeks. Performing additional monthly treatments, after completing the first 7 treatments, will help you maintain your smooth, silky, hair-free results.

### treatment durations

For effective results, mē should be used according to the following guidelines:

1 minute  
per underarm

2 minutes  
per bikini area

3 minutes  
per arm

10 minutes  
per full leg

2 minutes  
per shoulder

5 minutes  
per abdomen

6 minutes  
per chest


12 minutes  
per back

## quick reference guide


**Note:** Always start at LOW elōs level. If you feel comfortable, increase the elōs setting for higher efficacy

Contact the experts at 1-877-651-1177 to help get you started and set up your treatment routine. Our goal is to help you get the best results possible!

### done with mē?

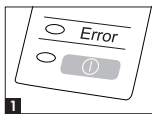1. Turn off the mē system by pressing the power button on the base unit.

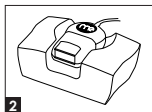2. Place the applicator into its cradle in the base unit and disconnect unit from power outlet.

3. Clean the cartridge light output window with a non-alcoholic based wet wipe or soft tissue.

4. Store the system in a cool and dry place AWAY from the reach of children.

5. Gently rinse your skin with lukewarm water and pat dry.

6. Gently apply moisturizing cream or aloe Vera gel to your skin.

### your mē package includes:

- mē system powered by elōs technology
- one elōs cartridge (already installed)
- accessory protective cap (already installed)
- elōs protective cap
- power cord
- user guide including a quick reference guide

## quick reference guide

### mē accessories (sold separately):

- replacement elōs cartridge
- pouch for storage

For more details about accessories, go to [www.mepower.com](http://www.mepower.com)

### device components (Figure 1):

1. base unit.
2. cradle.
3. control panel with:
  - power button
  - elōs energy setting buttons
  - indicators for the percentage (%) of pulses remaining
  - accessory activation button
4. applicator (hand piece with cord).
5. elōs cartridge with light output window and RF contact bars.

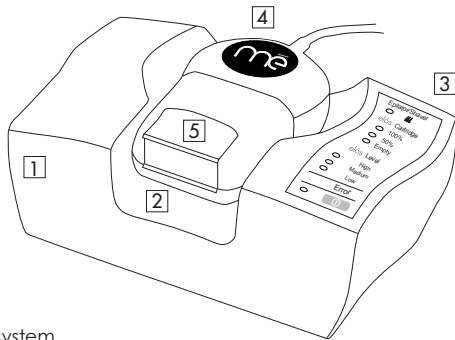
Figure 1: mē system

## quick reference guide

### applicator components (Figure 2):

6. light output window.

7. RF chrome contact bars (Both RF chrome bars must be in complete contact with the skin for the light pulses to flash).

8. elōs cartridge: comes attached to the applicator and is preloaded with elōs pulses. You will need to replace the cartridge when the Empty indicator light on the base unit is red. \*

9. elōs activation button.

10. cooling outlet.

11. accessory protective cap (not to be removed).

\*If you don't have a replacement cartridge – make sure you place the elōs protective cap instead of the used cartridge. This is in order to protect the inner parts of applicator.

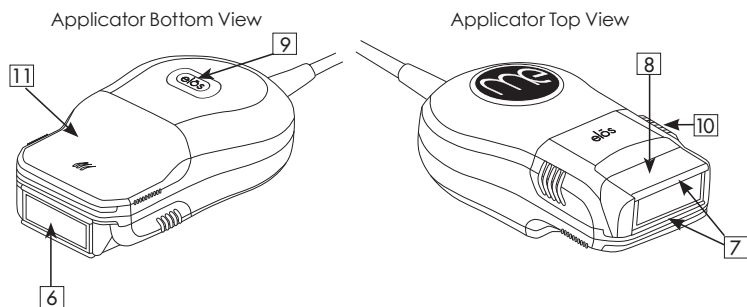
Figure 2: Top and bottom view of applicator with cartridge

## quick reference guide

### at-home professional hair removal

Every hair in our body goes through a complete growth cycle, from an active to resting phase. Light-based treatments work best on hairs in their active phase. How long each cycle lasts varies from person to person and between the different parts of the body. Other factors such as genetics, metabolism and age also influence the duration of each growth stage.

Hair reduction is a gradual and cumulative process. Our clinical studies have shown that there is less hair growth by using mē. How long it takes to see visible results and how long the results last vary from person to person. Periodic treatments are recommended to maintain these results. Clinical studies have shown that when hair does grow back between each of the treatments, there will be less of it.

### important safety information

Review the user guide including the warnings and safety information provided before using mē. The quick reference guide is only supplemental and not intended to replace your user guide.

### customer service contact

Phone: Toll free (NA): 1.877.651.1177

Hours of operation: Mon-Fri, 9am-5pm Eastern Standard Time

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© January 2013 All rights reserved. Syneron and mē are trademarks of Syneron Beauty Ltd.