
## High Speed Rectal Thermometer

User's Guide TH060

Thank you for purchasing this Advanced Solutions™ thermometer. We're proud to be your choice to help you keep your child healthy and safe. Please be sure to read all warnings, instructions and notations before using this product.

### Basics of Your Thermometer

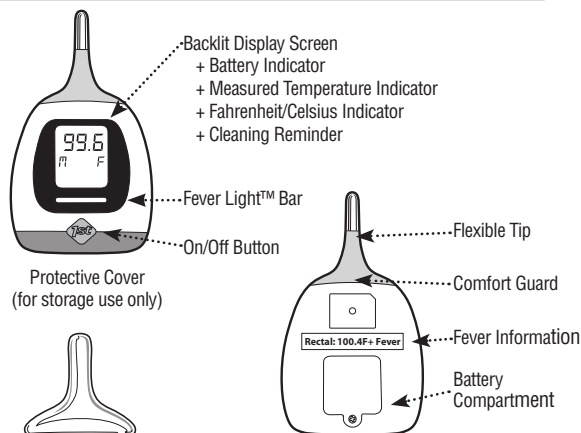Visit [www.safety1st.com](http://www.safety1st.com) and look for "Safety 1st TV" for information about how to take a rectal temperature.

### Fever Light™ Indicator

If your child has an elevated temperature (100.4°F +), the red FeverLight™ bar will light up and multiple beeps will sound. (See the Normal Body Temperature section for details about what constitutes a fever for your child.) High or prolonged fever requires medical attention, especially for young children. Contact your child's physician immediately.


### Normal Body Temperature

Temperature readings are different from person to person and also vary depending upon the body location from which the temperature is taken. A rectal temperature tends to be approximately 1° higher than an oral temperature for the same individual.

"Normal" temperature refers to each individual's average temperature while they are well. You may want to measure and record your baby's temperature a few times when he or she is well so that you will know what is normal for your own child.

On average, temperature ranges are as follows:

| TEMPERATURE-METHOD | NORMAL TEMPERATURE | FEVER RANGE |
|--------------------|--------------------|-------------|
| Rectal             | 98.0° to 100.3°F   | 100.4°F +   |
| Oral               | 97.0° to 99.5°F    | 99.6°F +    |
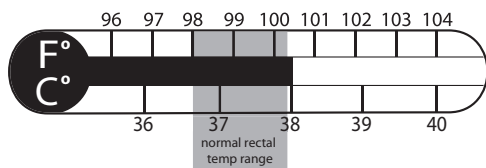
## WARNING

- This is **NOT** an ear thermometer. **NEVER** put probe into an ear.
- This is a rectal thermometer **ONLY**. It should **NEVER** be used for oral, underarm, ear or any other mode of temperature measurement.
- **NEVER** allow a child to play with this thermometer. It is not a toy. This is a sensitive electronic device and may be damaged if not used properly under parental supervision.
- **ALWAYS** supervise your child when using thermometer. **NEVER** allow your child to use it alone or without close parental supervision.
- Avoid sudden impact and dust.
- Do not expose the thermometer to extreme temperatures, humidity, direct sunlight or shock.
- Replace only with CR1632 or equivalent lithium battery. Use of other types of batteries may present a risk of fire or rupture.

## BATTERY CAUTION


- Keep batteries out of child's reach. Swallowing a battery could be fatal.
- Used batteries should be disposed of properly.
- Battery should not be charged or placed into extreme heat as it may explode.
- Remove batteries when stored for long periods of non-use.

### IMPORTANT INFORMATION

- Read these instructions carefully to ensure accurate temperature readings. Temperature readings are affected by many factors including exercise, drinking hot or cold beverages prior to taking a reading as well as technique.
- Rectal temperatures vary by approximately 0.9°F/0.5°C.
- **DO NOT** allow child to walk or run during temperature taking. Probe could present a hazard.

### How to Measure Temperature

*Please read this section carefully and all the way through before using the thermometer. Ensure that thermometer is at room temperature, waiting 15 minutes if necessary.*

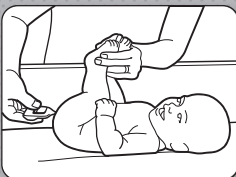
- Before first use, pull plastic tab to turn on thermometer. On/Off button to turn the power on. (One beep will sound to indicate it is on.)
- The last measurement taken will appear for 3 seconds.
- When the thermometer is ready to take a measurement, the bar will flash from left to right and "TAKE TEMP" will appear on the screen.
- To change measurements from Fahrenheit to Celsius, hold the On/Off button for 5 seconds. The screen will display "--.C" or "--.F" Release the button and press it again to select the correct mode. The thermometer will be ready to use after 3 seconds.
- Lubricate probe tip with a water-soluble lubricant or petroleum jelly.
- For easier insertion, your child should be positioned with legs at a 90-degree (L-shaped) angle – baby lying on back with legs in the air.

- Gently insert the probe tip no more than ½ inch (1.3 cm) into the rectum. Do not force the tip into the rectum if resistance is encountered. This thermometer's Comfort Guard will help prevent over-insertion.
- Once you have positioned the thermometer, hold it steady while the temperature is taken. The thermometer will make one long beep, the back light will illuminate and the "Clean Me!" icon will blink when the measurement is complete.
- If your child's temperature is at or above 100.4°F, the red Fever-Light™ bar will light up for 8 seconds and multiple beeps will sound, indicating that your child has a fever.
- Press the On/Off button to turn the thermometer off or it will shut off automatically after about 1 minute.
- Note: Wait at least 1 minute before taking another measurement.
- Clean the thermometer before putting it away. (See "Care and Maintenance")

### Care and Maintenance

- Clean the probe and tip with an alcohol swab before and after each measurement.
- Use a soft, dry cloth to clean thermometer body. Never use abrasive cleaning agents, thinners or benzene for cleaning.
- Do not scratch the surface of the display.
- Do not expose the thermometer to extreme temperatures, humidity, direct sunlight or shock.
- Store at room temperature.

### Normal Body Temperature (CR1632 or Equivalent Button Cell)

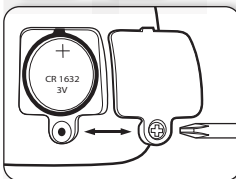To alert you when the battery should be replaced for this thermometer, the battery icon will display half of the battery and will flash. However, with half of the battery bar left, the thermometer can still be used for at least 20 more measurements. When the battery is completely empty, it will turn off automatically.

#### TO REPLACE THE BATTERY

- Use a small screwdriver to loosen the screw to the battery compartment. (See "Know Your Unit" section.)
- Slide the cover off, remove the old battery and insert a new CR1632 or equivalent button cell battery.
- Slide the battery cover on again and secure the screw.
- Note: Properly dispose of the batteries, keeping them away from small children and heat.

### If You See an "ERR" message

The thermometer screen will display "Err" if there is a system malfunction. Try turning it off and turning it on again. The thermometer will turn off automatically in about 1 minute if it is not manually turned off.

### Technical Specifications

**MEASUREMENT RANGE:** 90°F to 109°F

**MEASUREMENT ACCURACY:** 0.2°F between 96°F to 107°F

#### TEST PROTOCOL:

- Put the tip of thermometer into water bath, approximately 1/2" below the surface.
- Wait for 10 seconds, then turn on the thermometer by pressing the ON/OFF button. The thermometer is now in standard mode.
- When the measurement is completed, record the temperature reading.

### One (1) Year Limited Warranty

Dorel Juvenile Group, Inc. warrants to the original purchaser that this product (Thermometer), is free from material and workmanship defects when used under normal conditions for a period of one (1) year from the date of purchase. Should the product contain defects in material or workmanship Dorel Juvenile Group, Inc., will repair or replace the product, at our option, free of charge. Purchaser will be responsible for all costs associated with packaging and shipping the product to Dorel Juvenile Group Consumer Relations Department at the address noted on the front of this document and all other freight or insurance costs associated with the return. Dorel Juvenile Group will bear the cost of shipping the repaired or replaced product to the purchaser. Product should be returned in its original package accompanied by a proof of purchase, either a sales receipt or other proof that the product is within the warranty period. This warranty is void if the owner repairs or modifies the product or the product has been damaged as a result of misuse.

This warranty excludes any liability other than that expressly stated previously, including but not limited to, any incidental or consequential damages.

**SOME STATES DO NOT ALLOW THE EXCLUSION OR LIMITATION OF INCIDENTAL OR CONSEQUENTIAL DAMAGES, SO THE ABOVE LIMITATION OR EXCLUSION MAY NOT APPLY TO YOU.**

**THIS WARRANTY GIVES YOU SPECIFIC LEGAL RIGHTS, AND YOU MAY ALSO HAVE OTHER RIGHTS THAT MAY VARY FROM STATE TO STATE.**

©2013 Cosco Management, Inc.  
All Rights Reserved. Todos derechos reservados. [www.djgusa.com](http://www.djgusa.com)  
(800) 544-1108 [www.safety1st.com](http://www.safety1st.com)  
Made in CHINA. Hecho en CHINA.  
Styles and colors may vary. Los estilos y los colores pueden variar.  
Distributed by (distribuido por) Dorel Juvenile Group, Inc.,  
2525 State St., Columbus, IN 47201-7494  
Dorel Distribution Canada, 873 Hodge, St. Laurent, QC H4N 2B1  
01/25/13 4358-6051A

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