mē™

quick reference guide
guide de référence rapide
guía de referencia rápida

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thank you!

Congratulations on purchasing your mē system!

mē is an innovative device that can be used in the comfort and privacy of your own home to permanently reduce unwanted hair.

mē is the only at-home permanent hair reduction device to use patented elōs technology and is clinically proven to be safe and effective for all skin tones. Gentle, easy to use and just minutes a treatment, now you can experience the freedom of silky, smooth skin without the hassles of shaving and waxing.

Before using mē, visit our web site at www.mepower.com and make sure to register your mē and watch our demonstration video.

at-home professional hair removal

Every hair in our body goes through a complete growth cycle, from an active to resting phase. Light-based treatments work best on hairs in their active phase. How long each cycle lasts varies from person to person and between the different parts of the body. Other factors such as genetics, metabolism and age also influence the duration of each growth stage.

Hair reduction is a gradual and cumulative process. Our clinical studies have shown that there is less hair growth by using mē. How long it takes to see visible results and how long the results last vary from person to person. Periodic treatments are recommended to maintain these results.

Clinical studies have shown that when hair does grow back between each of the treatments, there will be less of it.

the science behind mē and elōs (Electro-Optical Synergy)

From Syneron and the scientists that invented IPL comes the new breakthrough patented hair removal technology called elōs. This unique combination of Intense Pulsed Light (IPL) and Radio Frequency (RF) is used to remove hair through selective targeting of the hair follicle in a safe yet effective manner. elōs is clinically proven to deliver unprecedented safety, speed and effectiveness for all skin tones.

the benefits of mē

- suitable for all skin tones
- suitable for all hair colors, even fair hair (epilator is recommended)
- for both body and face
- fast and easy, 10 minutes for full leg treatment and only 32 minutes for full body treatment
- very safe thanks to the effective elōs technology
- for both men and women
- visible results within 6-8 weeks
- dedicated face adaptor available, for safe and effective facial treatment
- no preparations required when using the dedicated epilator or shaver
- greatest number of pulses available in one cartridge!

important safety information

Review the user guide including the warnings and safety information provided before using mē. This quick reference guide is only supplemental and not intended to replace your user guide.

innovation and design

Syneron Beauty Ltd. Kochav Yokneam Bldg.
Yokneam Industrial Zone P.O Box 14 Yokneam Illit 20692, IL

customer service contact

contact us at www.mepower.com

innovation and design

Syneron Beauty Ltd. Kochav Yokneam Bldg.
Yokneam Industrial Zone P.O Box 14 Yokneam Illit 20692, IL

is mē for you?

mē is intended for use by women and men on body for facial hair below the neck and on ALL skin tones (skin color you are born with) and by women for facial use below the cheek line.

Body and face:

mē is intended for use on ALL hair colors.

Note: For better results, the epilator attachment is recommended for body treatment of naturally light blond, red and white/grey hair

Note: The dedicated mē accessory epilator should not be used on the face

your mē package includes:

- mē system powered by elōs technology
- one elōs cartridge (already installed)
- accessory protective cap (already installed)
- elōs protective cap
- power cord
- user guide including a quick reference guide
- protective eyewear

mē accessories:

Sold separately or included with your mē system

- shaver
- epilator
- face adaptor
- replacement elōs cartridge
- pouch for storage

For more details about accessories, go to www.mepower.com

device components (Figure 1):


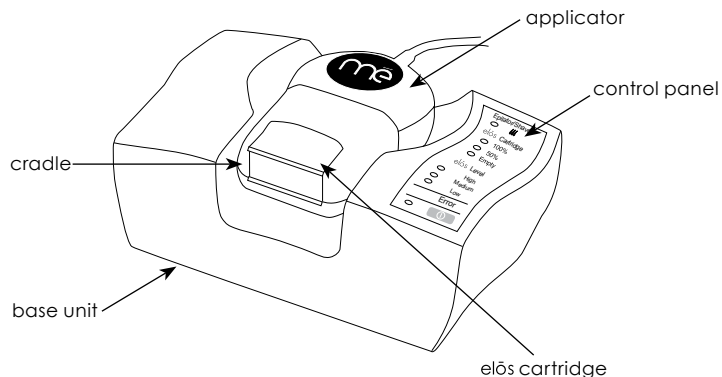
Figure 1: mē system

applicator components (Figure 2):

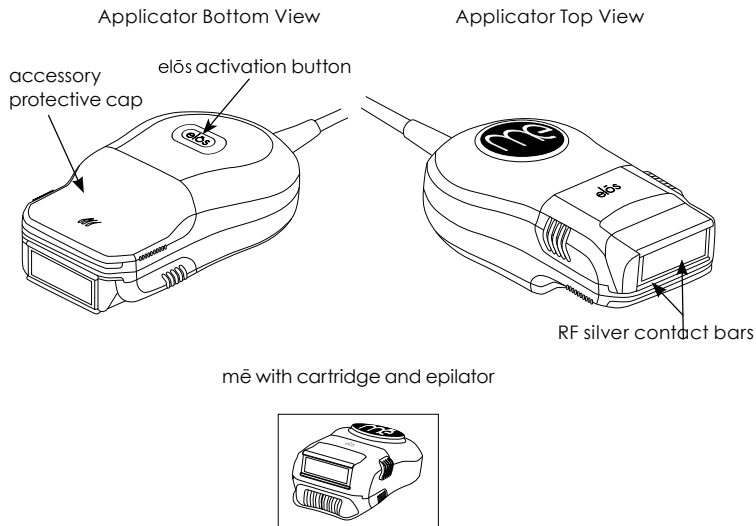*If you don't have a replacement cartridge – make sure you place the elōs protective cap instead of the used cartridge. This is in order to protect the inner parts of applicator.

Figure 2: Top and bottom view of applicator

how to use


step 1: preparing for mē

- Watch our demonstrational video at www.mepower.com
- Remove hair in the area you intend to treat with a razor or electric shaver, or alternatively, use the dedicated mē epilator/ shaver accessory


step 2: sensitivity test

Before treating any new area with mē, test a small patch of skin near the area you intend to treat.

step 3: treating with mē

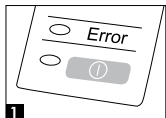
1

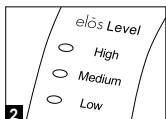2


3

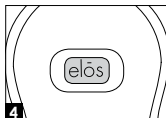4

Body:

1. Press the power button on the base unit to turn the system on.
2. Select the appropriate elōs level (Low- Medium or High). It is recommended to start on Low and if you feel comfortable use a higher elōs level.
3. If you are using the designated shaver/ epilator, connect it in place of the protective cap and press the accessory button on the control panel once to run the epilator mechanism at low speed. If desired, press it again to run the epilator at high speed. The shaver cartridge operates only at speed level 1.
4. Press the elōs activation button on the applicator.

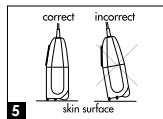5


5. Place the applicator on the area you want to treat, ensuring full contact of the two silver RF delivery bars with the skin. Do not press, just gently touch the skin. System will start emitting elōs pulses.

6. Glide the applicator slowly and continuously on the treatment area without stopping in one place.

Do not return to the same spot for at least 10 seconds.

Face:

Removing unwanted hair from the face is very different from removing unwanted hair from the body. Please read the user guide thoroughly before starting to remove facial hair using your mē. Be sure to watch the dedicated instructional video for face treatment on www.mepower.com.

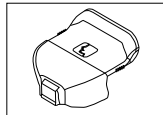When treating facial hair, it is recommended to use the **face adaptor**. The face adaptor has been designed to accommodate smaller and more delicate treatment areas on the face. It delivers energy in single pulse mode at the desired energy level.

The face adaptor may or may not be included with this package but can be purchased separately as an accessory. Face adaptor design may change.

Note: You may treat hair on the face from the cheek line and down


facial skin tones

mē is suitable for all skin tones. Make sure you perform a skin sensitivity test before you start your treatment (refer to the user guide for more info).


1

1. Clean your skin and remove the hair on the surface of the skin. You may use tweeze, wax, shave, etc. Do not use the mē epilator accessory on your face as it is not designed for facial use.

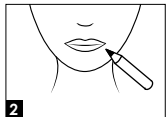2

2. Use white eyeliner to conceal areas that should not be treated such as the edges of your lips and facial moles.

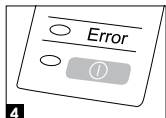4

3. Ensure you are in a well lighted area when performing your treatment.

4. Turn on your mē device and select the first elōs level on the base unit, which is Low. Put on the provided eye goggles or sunglasses.


4

With A Face Adaptor

5. Remove the accessory protective cap and connect the face adaptor. The face adaptor will allow the mē to emit only single pulses with the desired energy level.

6. Press the elōs cartridge activation button on the applicator to turn on the applicator. Gently position the applicator on the area you want to treat. Press once the elōs activation button to emit a single pulse.

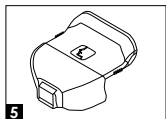
5

7. Reposition the applicator on the next area to be treated and press the elōs activation button again to emit another pulse. Repeat these previous steps until you have covered the entire facial area. Wait at least 10 seconds before you treat over the same spot.

Without A Face Adaptor

5. Press the elōs activation button on the applicator to turn on the applicator. Gently position the applicator on the area you want to treat. Upon contact - a pulse will be emitted. Remove the applicator after a single pulse.

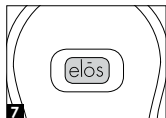7

6. Reposition the applicator on the next area to be treated, another pulse will be emitted, remove the applicator after a single pulse. Repeat these previous steps until you have covered the entire facial area. Wait at least 10 seconds before you treat over the same spot.

Note: Always start with LOW elōs level. If you feel comfortable, you may increase the elōs level for better results

Note: Without the face adaptor the system will emit continuous pulses, handling must be very cautious and following the instructions

Note: For men it is not advisable to treat the beard as it might create uneven results, which may produce changes in facial characteristics that may not be desirable

get your professionally smooth skin

For best results, use mē on the same area once a week for the first 7 times. Performing additional maintenance treatments, after the first 7 treatments, will help your hair removal results last longer.

* If you feel that the hair is too noticeable, you can treat that specific area every 2 days.

When treating the face, use mē twice a week for 2 weeks.

Then perform additional treatments every two weeks to help results last longer.

1st week 2nd week treatment durations

For effective results, should be used according to the treatment times indicated below.

The minimum treatment times recommended below may vary from person to person. We recommend to see that each treatment area is covered by 2-3 pulses:


2 minutes

1 minute
per underarm2 minutes
per bikini area3 minutes
per arm10 minutes
per full leg2 minutes
per shoulder5 minutes
per abdomen6 minutes
per chest12 minutes
per back

Note: Start your treatment at Low elös level.

If you feel comfortable, increase the elös setting for higher efficacy

Contact the experts at www.mepower.com to help get you started and set up your treatment routine. Our goal is to help you get the best results possible.

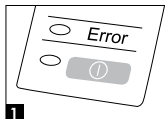
1

done with mē?

1. Turn off the mē system by pressing the power button on the base unit.

2. Place the applicator into its cradle in the base unit and disconnect unit from power outlet.

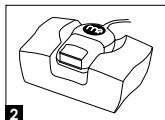2

3. Clean the cartridge light output window with a non-alcohol based wet wipe or soft tissue, clean accessories as well.

4. Store the system in a cool and dry place AWAY from the reach of children.

6

5. Gently rinse your treated areas with lukewarm water and pat dry.

6. You may gently apply moisturizing cream or aloe vera gel to your skin.

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© May 2013 All rights reserved. Syneron and mē are trademarks of Syneron Medical Ltd and may be registered in certain jurisdictions. Publication HU-PB00961.