beurer Aroma HK 65

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ENGLISH

Dear Customer,

We are happy that you have decided on a product from our range. Our name stands for high-quality and thoroughly tested products from the areas of warmth, blood pressure, body temperature, weight, gentle therapy, massage and air.

Please read these instructions for use carefully and follow the directions.

Best regards,
The Beurer Team

Items included in the package: 1 heating pad, 1 switch, 1 cover, 2 herb inserts, these instructions for use

1. Important safety information – read carefully and keep for later use

Failure to observe the information provided below can lead to personal injury or material damage (electric shock, skin burns, fire). The following safety and hazard information not only serves to protect your health and the health of others, it also protects the product from damage. Therefore, observe this safety information and pass on these instructions together with the product.

- Only use this heating pad for the purpose described in these instructions.
- This heating pad is not intended for use in hospitals.

- This heating pad is not intended for use by persons (including children) with restricted physical, sensory or mental capabilities and/or by persons lacking the required understanding of the equipment and how it is used, unless they are supervised by a person who is responsible for their safety.
- This heating pad should not be used with persons who are helpless, small children or persons insensitive to heat (e.g. diabetics, persons with skin changes related to illness or persons with scarred skin in the area where the pad is to be applied, after taking any pain-relieving medication or alcohol).
- Excessively long use can cause the skin to be burned.
- Caution! Never fall asleep while the heating pad is switched on.
- Do not apply the heating pad to parts of the body which are inflamed, injured or swollen. If in doubt, seek medical advice before use.
- Do not use the herb inserts if you have a known hypersensitivity or allergy to the herbs listed under point 2. If you have any doubts, seek advice from a medical professional before use.
- The electric and magnetic fields emitted by this heating pad can interfere with the operation of pacemakers. However, they are far below the permissible limits: electric field strength: 5,000 V/m max., magnetic field strength: 80 A/m max., magnetic flux density: 0.1 millitesla max. Therefore, please consult your doctor and the manufacturer of your pacemaker before using this heating pad.
- Heating pad
 - only connect to the voltage specified on the appliance,
 - only use in conjunction with switch types specified on the heating pad,
 - never operate without supervision,
 - never switch on when folded or compressed,
 - never wedge between objects,
 - never make any sharp creases,
 - never use with animals,
 - never use when damp.
- When in use, the electronic components in the switch of the heating pad make the switch slightly warm. The switch should therefore never be covered or placed on the heating pad when it is being operated.
- Do not pull, twist or make any sharp kinks in the cables.
- Do not stick pins or pointed objects into this heating pad.
- This appliance may not be used by children unless the switch has been preset by a parent or a supervisory person, or the child has been sufficiently instructed in how to use this heating pad safely.
- Children should be supervised to ensure that they do not use the heating pad as a toy.
- Please check this heating pad frequently for signs of wear and tear or damage. If there are signs indicating the above or if this heating pad was used improperly, it must be returned to the manufacturer or retailer before further use.
- If the mains connection cable of this heating pad is damaged, it must be replaced by the manufacturer, its customer service department or a similarly qualified person to avoid any risk.
- Repairs may only be carried out by specially trained personnel or at an authorised repair workshop, as special tools are required. Incorrect repairs may cause considerable danger to the user.
- Should you still have questions on using our appliances, please contact our customer service.

2. About aromatic applications

By warming aromatic medicinal herbs, the essential oils within are emitted into the air as vapours that can be inhaled.

2.1 Cold insert

Contents of the "cold" herb mixture: 1 g aniseed, 1 g thyme, 2 g peppermint.

The essential oils produced by these herbs are well-known to help clear the airways and reduce irritation of the throat.

2.2 Relaxation insert

Contents of the "relaxation" herb mixture: 2 g lavender, 1 g lemon balm, 1 g melissa

The essential oils produced by these herbs are well-known to have a relaxing effect.

3. Getting to know your instrument

3.1 Aloe vera cover

The removable textile cover of the heating pad has been treated with extract of aloe vera. As a result the textile cover is very kind to the skin and pleasantly soft.

The textile cover is machine-washable.

However, please note that the skin-friendly properties of the textile cover wash out after a number of wash cycles. A replacement textile cover can be ordered from the customer service department using item number 140.834.

3.2 Herb inserts

The heating pad comes complete with 2 different herb inserts. The relaxation insert and the cold insert.

Place one of the two inserts into the heating pad cover depending on the desired effect. Place the heating pad on the body, in a position that enables the essential vapours to be inhaled. If used daily, the inserts must be changed after approx. 1 month, because the essential oils will then have almost entirely evaporated. New inserts are available in shops.

The aromatic heating pad is suitable only for personal use, not for medical or commercial purposes.


4. Proper Use

This heating pad is designed for warming the human body. It allows you to apply targeted warmth. Heat promotes the circulation and relaxes the muscles. This heating pad is not intended for use in hospitals or for commercial applications. In particular, it may not be used on infants, small children, persons insensitive to heat, helpless persons or animals.

Our heating pads made of breathable microfleece are very adaptable and agreeable to the skin. They can also be used without their fabric cover. The detachable switch allows the pad to be cleaned in the washing machine (see cleaning instructions).

5. Appliance description

- 1 Heating pad
- 2 Plug-in coupling
- 3 Cover
- 4 Switch
- 5 Herb insert

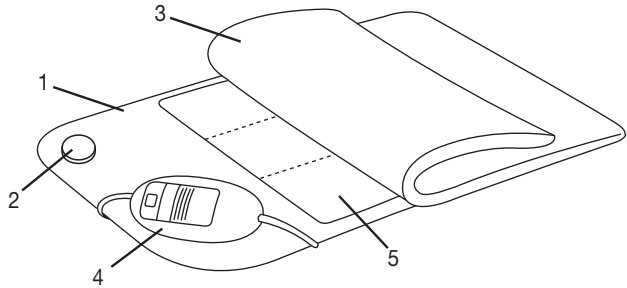6. Initial use

Remove the packaging from the desired herb insert and smooth it out onto the heating pad through the opened heating pad cover. Then close the zip fastening. Then check if the switch is connected to the heating pad via the plug-in coupling. Insert the mains plug into the socket. Place the heating pad on the upper body to enable the essential vapours to be inhaled.

7. Operation

7.1 Safety

This heating pad has a SAFETY SYSTEM. This electronic sensor technology prevents the entire surface of the heating pad from overheating by automatically switching off in the event of a fault. If the safety system switches off the heating pad automatically as a result of a fault, the indicator lamp is no longer illuminated even when the heating pad is switched on. Please note that, after a fault has occurred, the heating pad can no longer be used for safety reasons and must be sent to the service address indicated.

7.2 Selecting temperature

After switching on, the switch position will remain continuously illuminated. The fastest way to heat up the heating pad is to first set the switch to the highest setting. A lower setting can be selected later if required.

We recommend the lowest temperature setting for extended use of the heating pad.

7.3 Temperature settings

- Setting:
- 0 Off
 - 1 Minimum heat
 - 2 Medium heat
 - 3 Maximum heat

7.4 Automatic switch-off function

This heating pad is equipped with an automatic switch-off function. This turns off the heat supply approx. 90 minutes after the initial use. The switch position will then begin to flash. In order to reactivate the heating pad, the position switch must initially be set to the "0" position. After about 5 seconds it is possible to switch it on again.

If the heating pad will not be used after timed switch-off, it should be switched off (position "0") and the plug removed from the mains socket.

7.5 Fast heating

These heating pads have a fast heating function which allows them to heat up quickly within the first 10 minutes.

8. Cleaning and Care

Before cleaning the heating pad, always unplug it from the mains.

Smaller stains on the heating pad can be removed using a cloth or moist sponge together with a little liquid detergent for delicate fabrics if necessary.

The aromatic inserts must not be washed. Do not use any cleaning agents containing solvents on the appliance.

Please note that the heating pad should not be dry-cleaned, wrung out, machine-dried, pressed or ironed.


Before washing, first remove the plug-in connection together with the switch from the heating pad (see illustration). Set the washing machine to a gentle washing cycle at 30°C.

For ecological reasons, we recommend that the heating pad is only washed together with other textiles. Use a washing powder or liquid for delicate fabrics and dose it according to the manufacturers instructions.

Please note that excessive washing will wear out the heating pad. For this reason, the heating pad should only be machine-washed a maximum of 10 times throughout its entire life.

When drying the heating pad, do not use clothes pegs or anything similar, and never switch the heating pad on to dry it!

Do not reconnect the switch to the heating pad until the plug-in connection and the heating pad are completely dry.

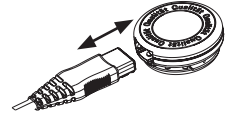9. Storage

If you are not going to use the heating pad for some time, we recommend storing it in its original packaging in a dry place without weighing it down.

Allow the heating pad to cool down beforehand.

10. Disposal

Follow local regulations regarding the disposal of materials. Please dispose of the blanket in accordance with the directive 2002/96/EG – WEEE (Waste Electrical and Electronic Equipment). If you have any queries, please refer to the local authorities responsible for waste disposal.

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